



EMSL ANALYTICAL, INC.

200 ROUTE 130 NORTH
CINNAMINSON, NJ 08077
PHONE: (845) 469-8671
FAX: (845) 231-6017

29 October 2013

Douglas Kent
TechLaw, Inc.
ESAT Region 8
16194 W. 45th Drive
Golden, CO 80403

Mr. Kent,


Thank you for your visit to our Libby laboratory on September 25th and 26th. Please find enclosed documentation of the corrective actions we discussed. These corrective actions include:

- 1) All EMSL/Libby analysts assign NaK codes to the 40 spectra of LA from the mine. Results for both before and after training are provided in a tabular format.
- 2) Index 3 diffraction patterns of NIST SRM 1867 actinolite. Include print-outs of EDS spectra for the 3 structures that had their diffraction patterns indexed.
- 3) Index 3 diffraction patterns of diopside from the bulk sample provided by ESAT. Include print-outs of EDS spectra for the 3 structures that had their diffraction patterns indexed.
- 4) The results of the intra-EMSL inter-lab study, which will include evaluation of NaK concordance.

Libby Amphibole EDXA interpretation

As you are aware LA particles are characterized by the presence of sodium, magnesium, potassium, calcium and iron. Aluminum is usually absent, but may occur at low levels in some structures. The use of EDS to classify an amphibole particle as LA is complicated by a number of factors, including a) inherent variability in chemical composition of different LA particles, b) dependence of the spectrum on random variables such as particle thickness, orientation, and proximity to other particles, c) variation between instruments in sensitivity to different elemental constituents, d) the presences of pyroxenes in the Rainey Creek Geologic Complex and e) subjective nature of spectra interpretation. For these reasons, it is not possible to define a single unique EDS spectrum for LA particles and an extensive library of references is needed to confirm a structure's identity. As such, the lab ran 40 spectra from the Libby Amphibole standards. These standards were used for the Libby Amphibole Spectra study in 2002. According to the SRC (Syracuse Research Center) these standards are from the following locations and contain the following minerals.

Sample Number	Source	Primary Amphibole Components
SW01201	Libby mine (site 20)	Tremolite, Actinolite
SW01281	Libby mine (site 28)	Richterite, Winchite
SW01231	Libby mine (site 23)	Tremolite, Actinolite, Richterite, Winchite
SW01TR1	NIST	Tremolite

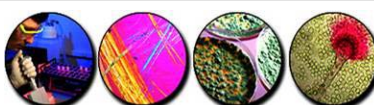

www.emsl.com



EMSL ANALYTICAL, INC.

200 ROUTE 130 NORTH
CINNAMINSON, NJ 08077

PHONE: (845) 469-8671
FAX: (845) 231-6017

Each analyst was instructed to independently assign mineral classification and NaK content for each spectra. I then personally reviewed with each individual analyst their data and spectra interpretations. It was

apparent from these discussions that the Libby staff had been trained to interpret that any energy channel above background should be considered clearly present. Due to the background and noise that are inherent to all EDXA systems this resulted in the majority of their spectra being classified as NaK. Additional training on EDXA theory and the production of statistically significant peaks was performed. After our discussion and the discussion with Doug Kent and Nathan DelHirro the exercise was run again. The results of these exercises are attached as well as the training record that was provided to each TEM analyst at Libby. **Please see attachments EMSL1.**

SAED/EDXA study of NIST 1867 Actinolite

The lab prepared a TEM drop mount of the NIST 1867 Actinolite for analysis. The TEM analysts were asked to provide 3 spectra and 3 indexed diffraction patterns from fibers they encountered during analysis. **Please see attachments EMSL2.**

SAED/EDXA study of Diopside (Herschel, Ontario)


The Diopside sample was obtained from Ward's Science and provided to EMSL Libby and Cinnaminson courtesy of Doug Kent. Each TEM analyst was asked to provide 3 spectra and 3 indexed patterns of the Diopside. **Please see attachments EMSL3.**

Inter Lab EMSL QC

Seven samples were picked for QC. The samples were analyzed in accordance with Lab mod 29. The results of the inter-lab EMSL QC exchange are summarized below

EMSL Sample Number	Customer sample ID	Reprep Result	Interlab Result	Result
271300789-0004	4R-02146	1 NAM	1 NAM	Concordant
271300800-0002	4R-02123	1 FP	1 NAM	Discordant
271300675-0010	P3-00774	1 NC LA	1 NC LA	Concordant
271300675-0003	P3-00880	ND	ND	Concordant
271300573-0001	SM-10080	ND	1 NC LA	Concordant
271300838-0001	4R-01595	ND	ND	Concordant
271300764-0002	4R-02012	ND	ND	Concordant

The majority of the samples were concordant- two samples showed documentation errors and in one of those instances lead to a false positive. One sample was improperly counted at corporate and is

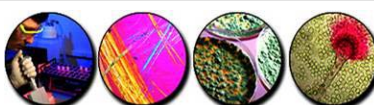

www.emsl.com



EMSL ANALYTICAL, INC.

200 ROUTE 130 NORTH
CINNAMINSON, NJ 08077

PHONE: (845) 469-8671
FAX: (845) 231-6017

discussed below.


For sample 271300573-0001 the incorrect fiber size recording rules were listed on the analytical worksheet. The analysis was for Low Magnification, PCME fibers only. However, the minimum fiber size listed on the analytical worksheet was 0.5 um.

EMSL04 analyzed the sample at high magnification (~ 20,000 xs) and found 1 structure. When compared to the re-prep analysis it was noted that the sample was analyzed at 4800x and it became apparent that the sizing information was incorrect on the benchsheets provided. This sample was part of the wildfire scenario which was analyzed at EPAs request prior to the SAP being finalized. The Libby Lab did receive written direction from Christina Progross on how to proceed with analysis. However these instructions were not forwarded to Corporate during the internal QC exchange. The SAP summary that was forwarded stated that the recording rules for the sample could be found in the "OU3 ISO modification #3: Rapid Air Analysis". This modification could not be found in either the Libby e-room or the OU3 e-room. This instance has been used as a training session with the staff to ensure that at receiving all the counting rules and pertinent information is reviewed before the samples get back to the lab and once back to the lab analysis should not commence until the analyst has reviewed the SAP summary and all the referenced documents in the SAP summary.

For Sample 271300800-0002 1 LA structure was listed on the bench sheet and Verified GO Map in GO C2: C8. The EDXA spectra provided documented that the structure was located in GO C1:G8. This structure was not located during the inter-lab verified analysis in either grid opening listed on the bench sheet/verified map or EDXA. During the reconciliation, the structure was located in GO C2:G8. The improper documentation of this grid opening led to a false positive being recorded. Roy Pescador investigated and determined that the root cause was not completing the analytical worksheet correctly prior to moving to the next grid opening. A corrective action was implemented to ensure that the correct grid openings are recorded on all documentation. Going forward the analyst will ensure that directly upon moving to a new grid, the grid ID will be documented before proceeding to the next grid opening. Prior to recording the grid opening the analyst will switch to low magnification and ensure the GO ID. Once GO ID is confirmed the analyst will document it on the analytical worksheet. EMSL's corrective action form is attached to this letter. **Please see attachments EMSL4.**

If you need any further clarification, please don't hesitate to contact me.

Robyn Denton | *Special Projects Manager*

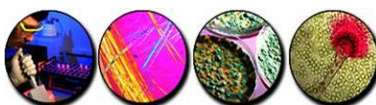


www.emsl.com



EMSL ANALYTICAL, INC.

200 ROUTE 130 NORTH
CINNAMINSON, NJ 08077
PHONE: (845) 469-8671
FAX: (845) 231-6017

ATTACHMENTS

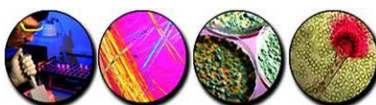

www.emsl.com




200 ROUTE 130 NORTH
CINNAMINSON, NJ 08077
PHONE: (845) 469-8671
FAX: (845) 231-6017

EMSL ANALYTICAL, INC.

Attachment EMSL1

September 25th, 2013 Pre-training

Source	Dominant Mineral Types	Spectra	R. Pescador				E.Wyatt-Pescador				K.Colberg				R. Denton			
			LA	NAM	Na	K	LA	NAM	Na	K	LA	NAM	Na	K	LA	NAM	Na	K
SW01201	Tremolite, Actinolite	B1	√		X	K	√		X	X	√		X	X	√		X	X
		B2	√		Na	X	√		Na	K	√		Na	X	√		X	X
		B3	√		X	K	√		X	X	√		Na	K	√		X	X
		B4	√		Na	K	√		Na	K	√		Na	X	√		X	X
		B5	√		X	K	√		X	K	√		Na	K	√		X	X
		B6	√		Na	X	√		Na	X	√		Na	K	√		X	X
		B7	√		X	K	√		Na	X	√		Na	X	√		X	X
		B8	√		Na	K	√		Na	K	√		Na	X	√		X	X
		B9	√		X	X	√		X	X	√		Na	X	√		X	X
		B10	√		Na	K	√		Na	K	√		Na	K	√		X	X
SW01231	Tremolite, Actinolite,Winchite, Richterite	D1	√		Na	K	√		Na	K	√		Na	K	√		Na	K
		D2	√		Na	K	√		Na	K	√		Na	K	√		Na	K
		D3	√		Na	K	√		Na	K	√		Na	K	√		Na	K
		D4	√		Na	K	√		Na	K	√		Na	K	√		Na	K
		D5		√	Na	K		√	Na	K		√	Na	K	√		Na	K
		D6	√		Na	K	√		Na	K	√		Na	K	√		Na	K
		D7	√		Na	K	√		Na	K	√		Na	K	√		Na	K
		D8	√		Na	K	√		Na	K	√		Na	K	√		Na	K
		D9	√		Na	K	√		Na	K	√		Na	K	√		Na	K
		D10	√		Na	K	√		Na	K	√		Na	K	√		Na	K
Boldface text denotes discordances																		

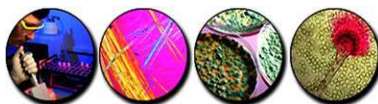

www.emsl.com



200 ROUTE 130 NORTH
CINNAMINSON, NJ 08077
PHONE: (845) 469-8671
FAX: (845) 231-6017


EMSL ANALYTICAL, INC.

Attachment EMSL1

September 25th, 2013 Pre-training

Source	Dominant Mineral Types	Spectra	R. Pescador				E. Wyatt-Pescador				K. Colberg				R. Denton			
			LA	NAM	Na	K	LA	NAM	Na	K	LA	NAM	Na	K	LA	NAM	Na	K
SW1281	Richterite, Winchite	E1	✓		Na	K	✓		Na	K	✓		Na	K	✓		Na	K
		E2	✓		Na	K	✓		Na	K	✓		Na	K	✓		Na	K
		E3	✓		Na	K	✓		Na	K	✓		Na	K	✓		Na	K
		E4	✓		Na	K	✓		Na	K	✓		Na	K	✓		Na	K
		E5	✓		Na	K	✓		Na	K	✓		Na	K	✓		Na	K
		E6	✓		Na	K	✓		Na	K	✓		Na	K	✓		Na	K
		E7	✓		Na	K	✓		Na	K	✓		Na	K	✓		Na	K
		E8	✓		Na	K	✓		Na	K	✓		Na	K	✓		Na	K
		E9	✓		Na	K	✓		Na	K	✓		Na	K	✓		Na	K
		E10	✓		Na	K	✓		Na	K	✓		Na	K	✓		Na	K
SW01TR1	NIST Tremolite	G1	✓		X	X	✓		X	X	✓		X	X	✓		X	X
		G2	✓		X	X	✓		X	X	✓		X	X	✓		X	X
		G3	✓		X	X	✓		X	X	✓		X	X	✓		X	X
		G4	✓		Na	X	✓		Na	X	✓		X	X	✓		X	X
		G5	✓		Na	X	✓		Na	X	✓		X	X	✓		X	X
		G6	✓		X	X	✓		X	X	✓		X	X	✓		X	X
		G7	✓		Na	X	✓		Na	X	✓		Na	X	✓		X	X
		G8	✓		X	X	✓		X	X	✓		X	X	✓		X	X
		G9	✓		X	X	✓		X	X	✓		Na	X	✓		X	X
		G10	✓		X	X	✓		X	X	✓		Na	X	✓		X	X

Boldface text denotes discordances

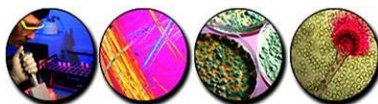

www.emsl.com



200 ROUTE 130 NORTH
CINNAMINSON, NJ 08077
PHONE: (845) 469-8671
FAX: (845) 231-6017


EMSL ANALYTICAL, INC.

Attachment EMSL1

September 26th, 2013-Post Training

Source	Minerals	Spectra	R. Pescador				E.Wyatt-Pescador				K.Colberg				R. Denton			
			LA	NAM	Na	K	LA	NAM	Na	K	LA	NAM	Na	K	LA	NAM	Na	K
SW01201	Tremolite, Actinolite	B1	√		X	X	√		X	X	√		X	X	√		X	X
		B2	√		Na	X	√		Na	X	√		Na	X	√		Na	X
		B3	√		X	X	√		X	X	√		X	X	√		X	X
		B4	√		X	X	√		X	X	√		X	X	√		X	X
		B5	√		X	X	√		X	X	√		X	X	√		X	X
		B6	√		X	X	√		X	X	√		Na	X	√		X	X
		B7	√		X	X	√		X	X	√		Na	X	√		X	X
		B8	√		X	X	√		X	X	√		X	X	√		X	X
		B9	√		X	X	√		X	X	√		X	X	√		X	X
		B10	√		X	X	√		X	X	√		Na	X	√		X	X
SW01231	Tremolite, Actinolite,Winchite, Richterite	D1	√		Na	K	√		Na	K	√		Na	K	√		Na	K
		D2	√		Na	K	√		Na	K	√		Na	K	√		Na	K
		D3	√		Na	K	√		Na	K	√		Na	K	√		X	K
		D4	√		Na	K	√		Na	K	√		Na	K	√		Na	K
		D5		√	Na	K		√	Na	K		√	Na	K	√		Na	K
		D6	√		Na	K	√		Na	K	√		Na	K	√		Na	K
		D7	√		Na	X	√		Na	X	√		Na	X	√		Na	X
		D8	√		Na	K	√		Na	K	√		Na	K	√		Na	K
		D9	√		Na	K	√		Na	K	√		Na	K	√		Na	K
		D10	√		Na	K	√		Na	K	√		Na	K	√		X	K

Boldface Text denotes discordances

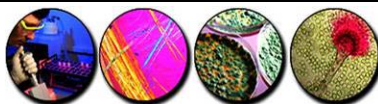

www.emsl.com



200 ROUTE 130 NORTH
CINNAMINSON, NJ 08077
PHONE: (845) 469-8671
FAX: (845) 231-6017


EMSL ANALYTICAL, INC.

Attachment EMSL1

September 26th, 2013-Post Training

Source	Minerals	Spectra	R. Pescador				E.Wyatt-Pescador				K.Colberg				R. Denton			
			LA	NAM	Na	K	LA	NAM	Na	K	LA	NAM	Na	K	LA	NAM	Na	K
SW1281	Richterite, Winchite	E1	√		Na	K	√		Na	K	√		Na	K	√		Na	K
		E2	√		Na	K	√		Na	K	√		Na	K	√		Na	K
		E3	√		Na	K	√		Na	K	√		Na	K	√		Na	K
		E4	√		Na	K	√		Na	K	√		Na	K	√		Na	K
		E5	√		Na	K	√		Na	K	√		Na	K	√		Na	K
		E6	√		Na	K	√		Na	K	√		Na	K	√		Na	K
		E7	√		Na	K	√		Na	K	√		Na	K	√		Na	K
		E8	√		Na	K	√		Na	K	√		Na	K	√		Na	K
		E9	√		Na	K	√		Na	K	√		Na	K	√		Na	K
		E10	√		Na	K	√		Na	K	√		Na	K	√		Na	K
SW01TR1	NIST Tremolite	G1	√		X	X	√		X	X	√		X	X	√		X	X
		G2	√		X	X	√		X	X	√		X	X	√		X	X
		G3	√		X	X	√		X	X	√		X	X	√		X	X
		G4	√		X	X	√		X	X	√		X	X	√		X	X
		G5	√		X	X	√		X	X	√		X	X	√		X	X
		G6	√		X	X	√		X	X	√		X	X	√		X	X
		G7	√		X	X	√		X	X	√		X	X	√		X	X
		G8	√		X	X	√		X	X	√		X	X	√		X	X
		G9	√		X	X	√		X	X	√		X	X	√		X	X
		G10	√		X	X	√		X	X	√		X	X	√		X	X

Boldface Text denotes discordances

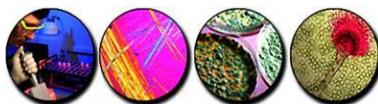




www.emsl.com




EMSL ANALYTICAL, INC.

200 ROUTE 130 NORTH
CINNAMINSON, NJ 08077
PHONE: (845) 469-8671
FAX: (845) 231-6017

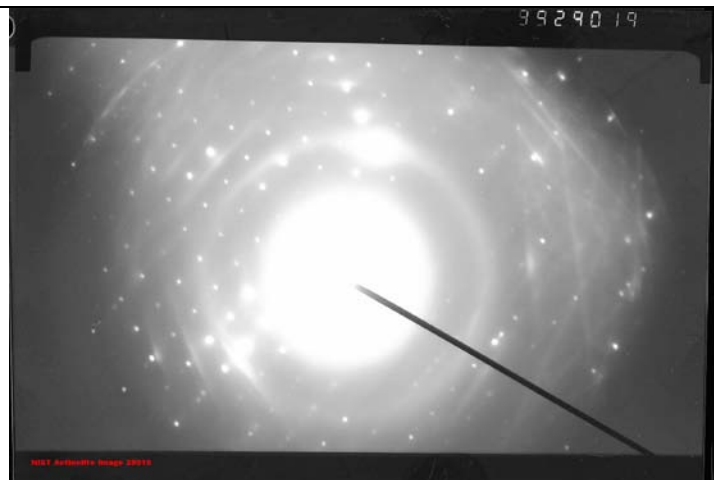Attachment EMSL2

	<p><i>SAED IMAGE 29019</i></p> <p><i>NIST Actinolite</i></p> <table><tr><td>Interrow</td><td>5.07 Angstroms</td></tr><tr><td>Hko</td><td>5.02 Angstroms</td></tr><tr><td>Hkl</td><td>4.76 Angstroms</td></tr><tr><td>Slant Vector Angle</td><td>68°</td></tr><tr><td>Zone Axis</td><td>[3 1 2]</td></tr></table> <p><i>Camera Constant 22.82 mm/Angstroms</i></p> <p><i>SAED Reference JCPDS</i></p>	Interrow	5.07 Angstroms	Hko	5.02 Angstroms	Hkl	4.76 Angstroms	Slant Vector Angle	68°	Zone Axis	[3 1 2]
Interrow	5.07 Angstroms										
Hko	5.02 Angstroms										
Hkl	4.76 Angstroms										
Slant Vector Angle	68°										
Zone Axis	[3 1 2]										
	<p><i>SAED IMAGE 29031</i></p> <p><i>NIST Actinolite</i></p> <table><tr><td>Interrow</td><td>5.37 Angstroms</td></tr><tr><td>Hko</td><td>9.13 Angstroms</td></tr><tr><td>Hkl</td><td>4.56 Angstroms</td></tr><tr><td>Slant Vector Angle</td><td>60°</td></tr><tr><td>Zone Axis</td><td>[1 0 0]</td></tr></table> <p><i>Camera Constant 22.82 mm/Angstroms</i></p> <p><i>SAED Reference JCPDS</i></p>	Interrow	5.37 Angstroms	Hko	9.13 Angstroms	Hkl	4.56 Angstroms	Slant Vector Angle	60°	Zone Axis	[1 0 0]
Interrow	5.37 Angstroms										
Hko	9.13 Angstroms										
Hkl	4.56 Angstroms										
Slant Vector Angle	60°										
Zone Axis	[1 0 0]										
	<p><i>SAED IMAGE 29035</i></p> <p><i>NIST Actinolite</i></p> <table><tr><td>Interrow</td><td>5.22 Angstroms</td></tr><tr><td>Hko</td><td>5.02 Angstroms</td></tr><tr><td>Hkl</td><td>2.56 Angstroms</td></tr><tr><td>Slant Vector Angle</td><td>85°</td></tr><tr><td>Zone Axis</td><td>[3 -1 -10]</td></tr></table> <p><i>Camera Constant 22.82 mm/Angstroms</i></p> <p><i>SAED Reference JCPDS</i></p>	Interrow	5.22 Angstroms	Hko	5.02 Angstroms	Hkl	2.56 Angstroms	Slant Vector Angle	85°	Zone Axis	[3 -1 -10]
Interrow	5.22 Angstroms										
Hko	5.02 Angstroms										
Hkl	2.56 Angstroms										
Slant Vector Angle	85°										
Zone Axis	[3 -1 -10]										

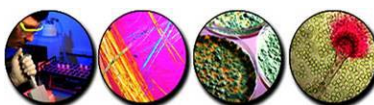


www.emsl.com



200 ROUTE 130 NORTH
CINNAMINSON, NJ 08077
PHONE: (845) 469-8671
FAX: (845) 231-6017


EMSL ANALYTICAL, INC.

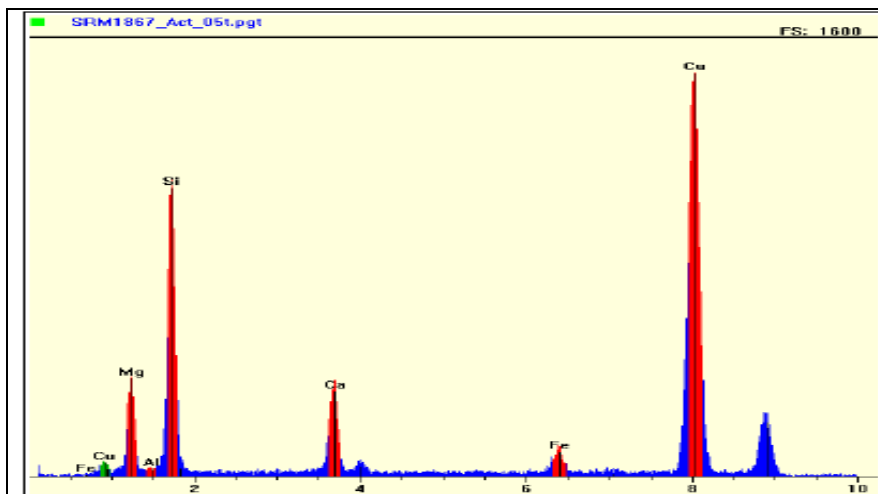

*EDXA Spectra NIST Actinolite
AC/XX*

Full scale 1600


*Corresponds to diffraction image
29019*

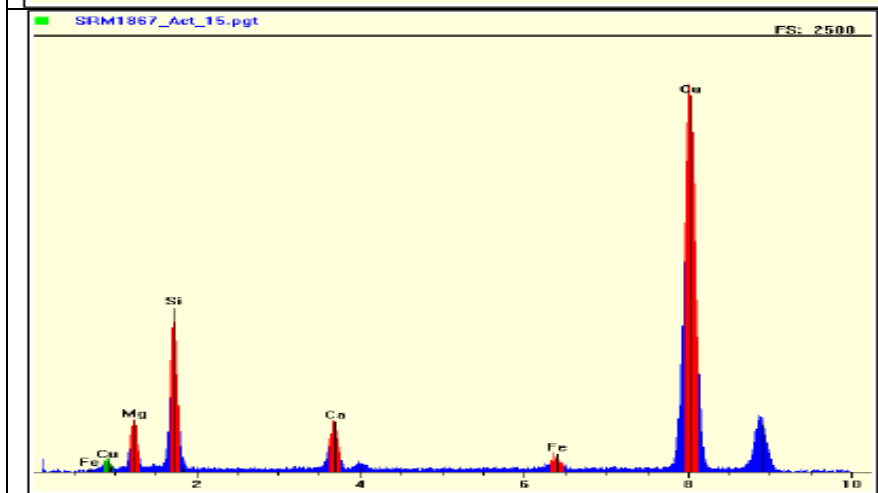

*EDXA Spectra NIST Actinolite
AC/XX*

Full scale 2500


*Corresponds to diffraction image
29031*

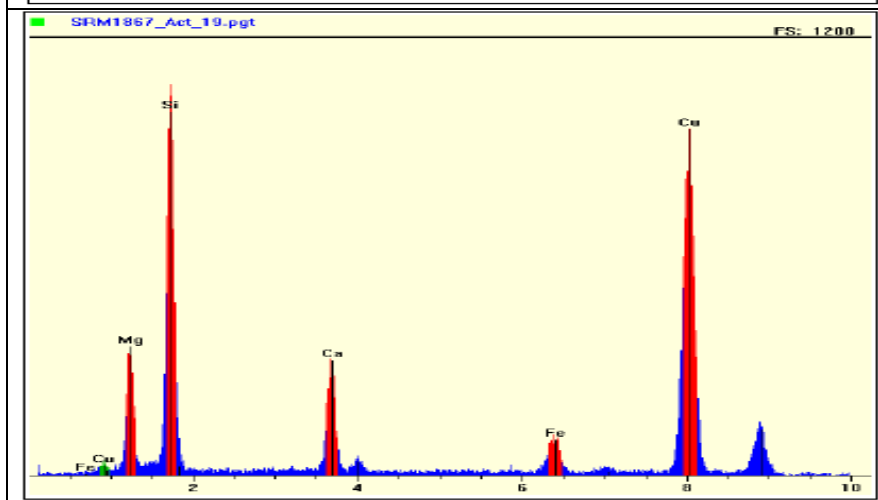

*EDXA Spectra NIST Actinolite
AC/XX*

Full scale 1200

*Corresponds to diffraction image
29035*

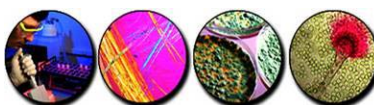




www.emsl.com




EMSL ANALYTICAL, INC.

200 ROUTE 130 NORTH
CINNAMINSON, NJ 08077
PHONE: (845) 469-8671
FAX: (845) 231-6017

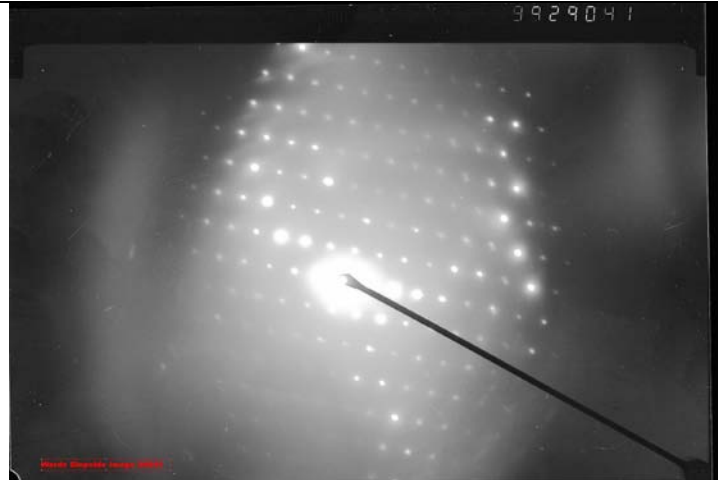Attachment EMSL3

	<p><i>SAED IMAGE 29041</i></p> <p><i>Wards Diopside</i></p> <table><tr><td>Interrow</td><td>5.19 Angstroms</td></tr><tr><td>Hko</td><td>6.29 Angstroms</td></tr><tr><td>Hkl</td><td>4.99 Angstroms</td></tr><tr><td>Slant Vector Angle</td><td>79°</td></tr><tr><td>Zone Axis</td><td>[-1 -1 0]</td></tr></table> <p><i>Camera Constant 22.82 mm/Angstroms</i></p> <p><i>SAED Reference JCPDS</i></p>	Interrow	5.19 Angstroms	Hko	6.29 Angstroms	Hkl	4.99 Angstroms	Slant Vector Angle	79°	Zone Axis	[-1 -1 0]
Interrow	5.19 Angstroms										
Hko	6.29 Angstroms										
Hkl	4.99 Angstroms										
Slant Vector Angle	79°										
Zone Axis	[-1 -1 0]										
	<p><i>SAED IMAGE 29051</i></p> <p><i>Wards Diopside</i></p> <table><tr><td>Interrow</td><td>5.27 Angstroms</td></tr><tr><td>Hko</td><td>3.68 Angstroms</td></tr><tr><td>Hkl</td><td>2.54 Angstroms</td></tr><tr><td>Slant Vector Angle</td><td>85°</td></tr><tr><td>Zone Axis</td><td>[-4 2 2]</td></tr></table> <p><i>Camera Constant 22.82 mm/Angstroms</i></p> <p><i>SAED Reference JCPDS</i></p>	Interrow	5.27 Angstroms	Hko	3.68 Angstroms	Hkl	2.54 Angstroms	Slant Vector Angle	85°	Zone Axis	[-4 2 2]
Interrow	5.27 Angstroms										
Hko	3.68 Angstroms										
Hkl	2.54 Angstroms										
Slant Vector Angle	85°										
Zone Axis	[-4 2 2]										
	<p><i>SAED IMAGE 29057</i></p> <p><i>Wards Diopside</i></p> <table><tr><td>Interrow</td><td>5.22 Angstroms</td></tr><tr><td>Hko</td><td>6.27 Angstroms</td></tr><tr><td>Hkl</td><td>2.56 Angstroms</td></tr><tr><td>Slant Vector Angle</td><td>80°</td></tr><tr><td>Zone Axis</td><td>[2 2 2]</td></tr></table> <p><i>Camera Constant 22.82 mm/Angstroms</i></p> <p><i>SAED Reference JCPDS</i></p>	Interrow	5.22 Angstroms	Hko	6.27 Angstroms	Hkl	2.56 Angstroms	Slant Vector Angle	80°	Zone Axis	[2 2 2]
Interrow	5.22 Angstroms										
Hko	6.27 Angstroms										
Hkl	2.56 Angstroms										
Slant Vector Angle	80°										
Zone Axis	[2 2 2]										

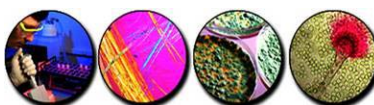


www.emsl.com



200 ROUTE 130 NORTH
CINNAMINSON, NJ 08077
PHONE: (845) 469-8671
FAX: (845) 231-6017


EMSL ANALYTICAL, INC.

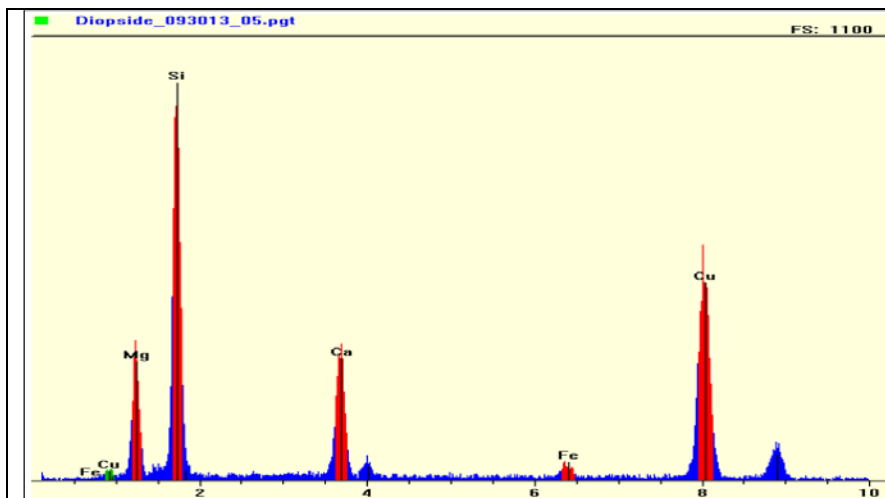

*EDXA Spectra Diopside Wards
PY/XX*

Full scale 1100


*Corresponds to diffraction image
29041*

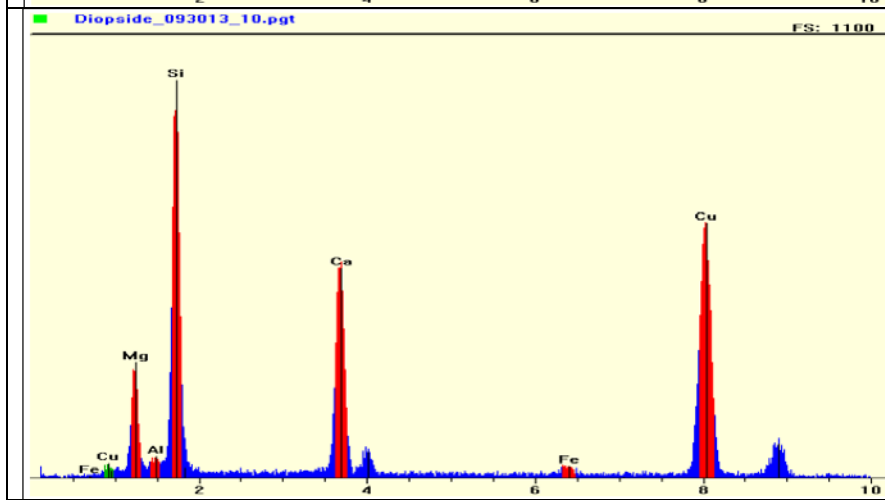

*EDXA Spectra Diopside Wards
PY/XX*

Full scale 1100


*Corresponds to diffraction image
29047*

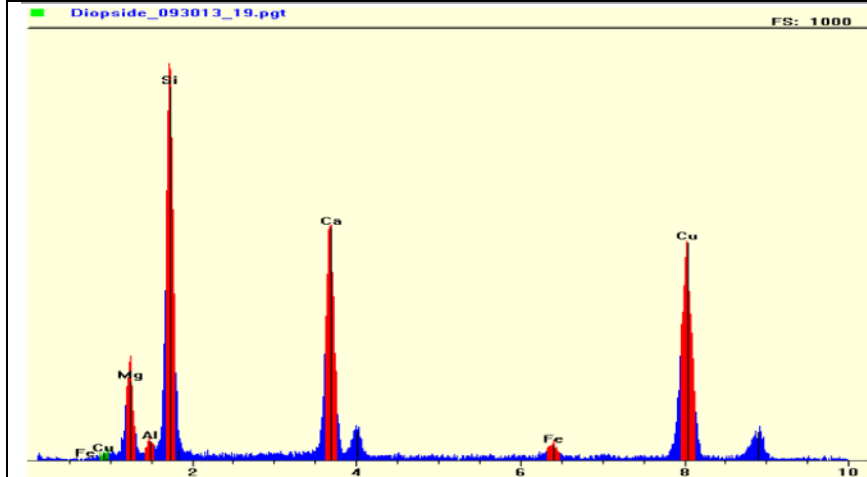

*EDXA Spectra Diopside Wards
PY/XX*

Full scale 1000

*Corresponds to diffraction image
29057*

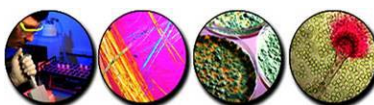

www.emsl.com




EMSL ANALYTICAL, INC.

200 ROUTE 130 NORTH
CINNAMINSON, NJ 08077
PHONE: (845) 469-8671
FAX: (845) 231-6017

SUPPORTING DOCUMENTATION

RAW EDXA SPECTRA NAK
LA SPECTRA INTERPRETATION TRAINING RECORDS
EDXA NIST ACTINOLITE
EDXA WARDS DIOPSIDE
SAED INDEXING FORMS
INTRA LAB QC RESULTS

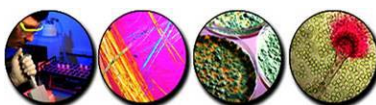

www.emsl.com



EMSL TRAINING RECORD

Employee: EXAMPLE

Title / Position:

Laboratory: Libby

Date of Training:

Trainer: Robyn Denton

Signature of Trainer:

Department: Asbestos

Topic: Libby Amphibole EDXA
Collection and
Interpretation

Details Of Training:

Robyn Denton visited the lab to discuss the collection of EDXA collection and Interpretation. During the training the Libby Amphibole standards were used. According to the SRC (Syracuse Research Center) these standards are from the following locations and contain the following minerals. These standards were used for the Libby Amphibole Spectra study in 2002.

Sample Number	Source	Primary Amphibole Components
SW01201	Libby mine (site 20)	Tremolite, Actinolite
SW01281	Libby mine (site 28)	Richterite, Winchite
SW01231	Libby mine (site 23)	Tremolite, Actinolite, Richterite, Winchite
SW01TR1	NIST	Tremolite

To Obtain a spectra:

Once the sample has been inserted, choose an intact grid opening.

- Attain eucentric position.
- Tilt to +32° to +35°
- Switch to an appropriate spot size. The appropriate spot size decision is made by taking into account the following.
 - The spot size should be able to produce a sufficient amount of x-ray counts to allow for adequate peak acquisition while not "flooding" the detector.
 - Sufficient x-ray counts are optimally 2000 counts per second. Change to a larger spot size if the counts are substantially lower than 2000 c/s.
 - "Flooding" the detector occurs when the amount of x-rays overwhelms the detector and the processor software producing a high dead time.
 - A high dead time is defined differently on different scope / detector combinations. In general:
 - For older scope / detector combinations "high" is when the deadtime exceeds about 20%
 - High deadtimes cause the detector to ignore a large percentage of x-rays that hit the window and the spectrum may be distorted.

Spectrum Collection Times:

- Collection time may vary depending on particle thickness, size and location. EMSL typically tries to collect until 1000 counts are reached in the Silicon peak.



- Note: The FS (full scale) count in PGT is referring to the highest peak in the spectrum, which when you are close to a grid bar could be Cu. So be sure to scale your spectrum appropriately.
- The longer the collection time the lower the background.
- Collection times <1000 counts may be used if the spectrum that is being generated has a low, flat background and distinct peaks in the Na and K regions.

To decrease the amount of deadtime, you can:

- Switch to a smaller spot size if the deadtime is substantially above the optimal values above.
- For PGT users, adjust the LLD . Under Xray setup, ADC tab, adjust the LLD to a value near 1.50. Collect a spectrum and see if the dead time has been lessened. Once you are comfortable with the settings, insert NIST Crocidolite and collect a spectrum. Determine your Na sensitivity.

Selecting a Spot Size:

- The size of the fiber/ particle will affect the choice of spot size. In general larger fibers (usually amphibole) will require the selection of a smaller spot size, while smaller fibers (chrysotile fibrils) require larger spot sizes.
- A trade off exists when choosing larger spot sizes.
 - The larger the spot size, the larger the area you are sampling.
 - If there are other particles present close to the particle of interest and you choose a large spot size, then the x-ray spectrum will exhibit characteristic x-rays from both the particle of interest and also nearby particles.

Interpreting/ Analyzing the spectrum

- The spectrum is a histogram, with the x-axis being the energy (KeV) and the Y-axis being the number of x-rays counted.
 - Each elemental x-ray will have a characteristic energy, that energy is measured by the detector.
 - The x-ray's energy is then plotted on the graph as one count, at the appropriate energy.
 - By plotting all x-rays, the graph will reveal peaks, energies that have a high numbers of counts, that can be identified as belonging to particular elements.
- Qualitative analysis consists of identifying the peaks present, using the software accompanying the detector processor. In general software can identify peaks in two ways:
 - By producing a line of interest and offering possible elemental matches for the line. The line can be moved to a peak of interest, possible elemental identifications for the peak are then offered.
 - Expected elements can be set, and their peak locations highlighted, so that any peak that happens to be present in this area is possibly from this element.
- Semi-quantitative analysis is the process of comparing these spectra, against known standard spectrum, to generate a possible matched identification.





Libby Amphibole Identification


LA particles are characterized by the presence of sodium, magnesium, potassium, calcium and iron. Aluminum is usually absent, but may occur at low levels in some structures. The use of EDS to classify an amphibole particle as LA is complicated by a number of factors, including a) inherent variability in chemical composition of different LA particles, b) dependence of the spectrum on random variables such as particle thickness, orientation, and proximity to other particles (back scatter), and c) variation between instruments in sensitivity to different elemental constituents, d) the presences of pyroxenes in the Rainey Creek Geologic Complex. Pyroxenes are chemically similar to LA qualitatively. Analysis of the Atomic formula can aid in providing identification, but require much data manipulation is only available on certain EDXA systems.

For these reasons, it is not possible to define a unique EDS spectrum for LA particles.


As noted above, there is substantial variability between particles and between instruments, but according to CHARACTERISTIC EDS SPECTRA FOR LIBBY-TYPE AMPHIBOLES dated March 15, 2005 and prepared by SRC the following peak magnitudes (relative to silicone) are typical for Libby Amphibole. These values should not be interpreted as absolute values for classifying LA; there is too much variability amongst particles and scopes. This should cover the majority of the LA encountered during analysis. Any structure that falls outside these bounds may still be LA, but will need more corroboration before being fully identified. (Morphology, diffraction) A spectrum from each standard is below. These spectra were collected on scope 27-2 in 2002.

Element	Typical Peak Magnitude		
	Mean	5th	95th
Na	0.36	0.15	0.64
Mg	3.08	2.47	3.49
Al	0.08	0.00	0.27
Si	10.00	10.00	10.00
K	0.25	0.06	0.53
Ca	1.81	0.96	2.63
Fe	1.08	0.57	1.88

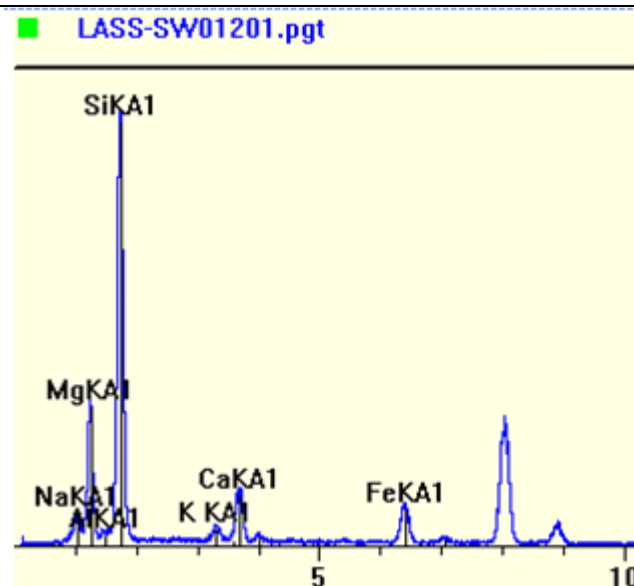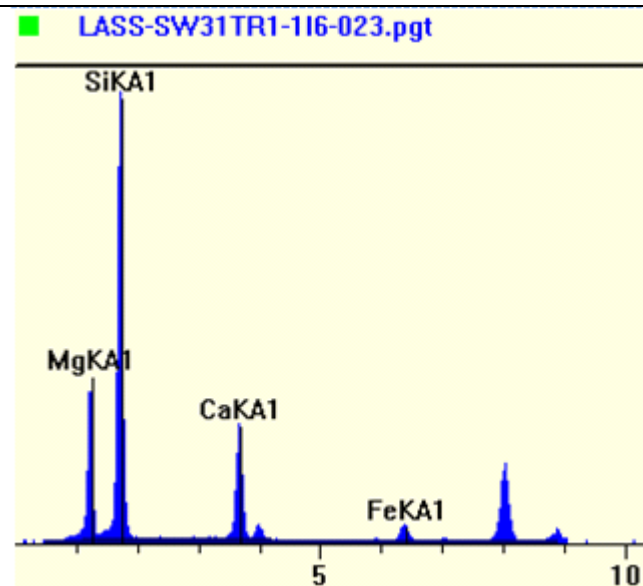

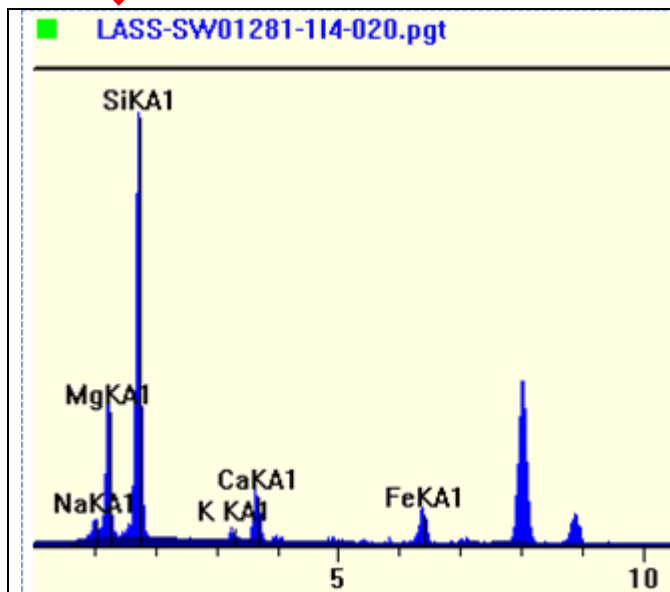

Primary Amphibole Components: Richterite, Winchite

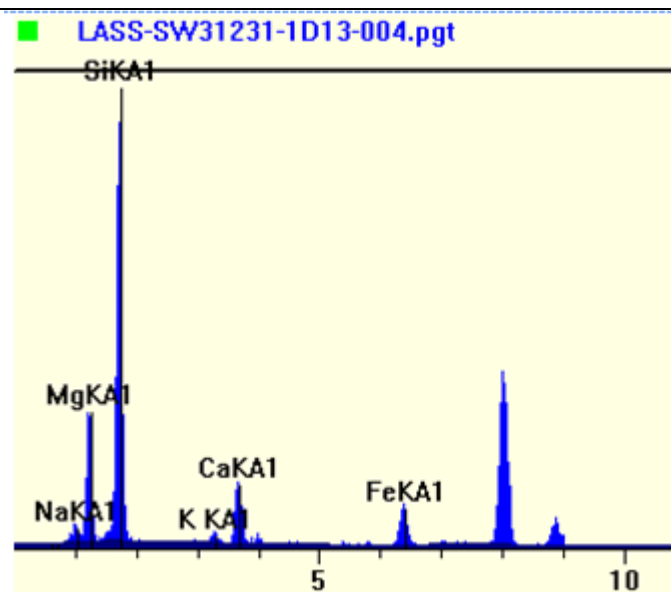


Primary Amphibole Components: Tremolite, Actinolite, Richterite, Winchite

To determine Na / K presence:

Per lab mod 66 C every structure identified as LA needs to have the presence of Na or K documented. In order to determine, an analyst will evaluate the spectra and determined if the Na or K is present in a statistically stable amount above background. A small peak should be present, thus reflecting that the element has a high number of counts and can be statistically valid. A ratio of a particular element to Silicon can be used and compared to the above the table. Lastly, care must be taken when observing Na, as the Cu L alpha overlaps the Na K alpha peak. We need to be certain that any Quantitative analysis exclude Cu L alpha counts.

Lab mod 66 E states that the Na or K must be clearly present. A little blip of K or Na is not clearly present and therefore should not be labeled as such.

Lab Mod 66 E has many spectra present that can be used as references. It is important to remember that these spectra were taken on different detector than what we use and our own detectors may vary slightly from the examples.

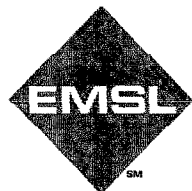

Energy Dispersive X-Ray Analysis

Qualitative Spectrum

EMSL ANALYTICAL, INC.

File: F:\Documen...ords\PGT Files\EMSL27-2\EMSL27-2 2013\sw01201\B01.pgt
Collected: September 24, 2013 08:31:54


Live Time: 75.62
Beam Voltage: 20.00


Count Rate: 2251
Beam Current: 2.00

Dead Time: 24.18 %
Takeoff Angle: 60.98

■ B01.pgt

FS: 1200

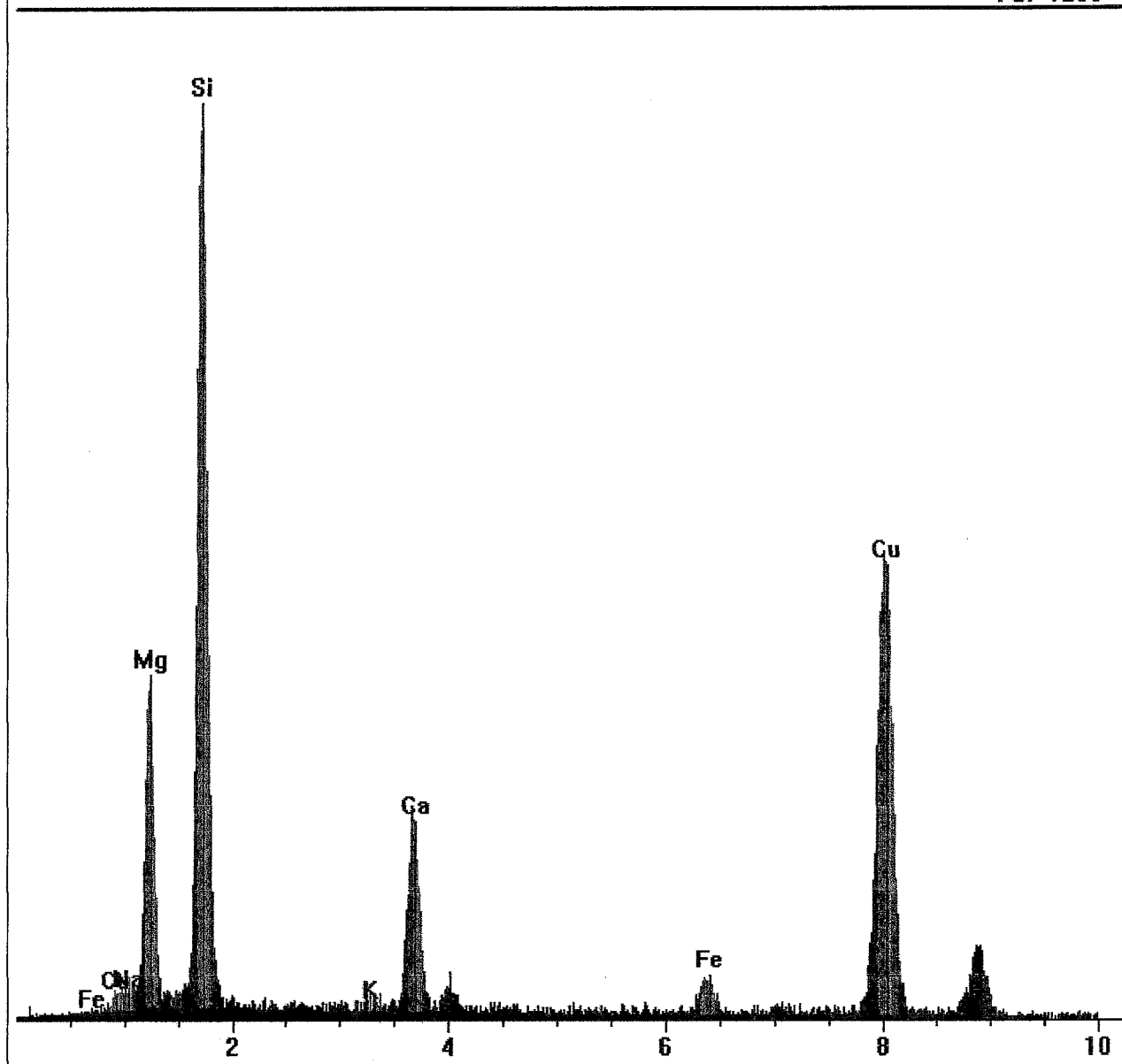

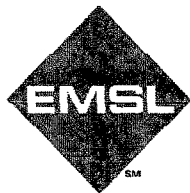

Energy Dispersive X-Ray Analysis Qualitative Spectrum

EMSL ANALYTICAL, INC.

File: F:\Documen...cords\PGT Files\EMSL27-2\EMSL27-2 2013\sw01201\B02.pgt

Collected: September 24, 2013 08:31:54

Live Time: 74.96

Count Rate: 1989

Dead Time: 21.99 %


Beam Voltage: 20.00


Beam Current: 2.00

Takeoff Angle: 60.98

■ B02.pgt

FS: 1000

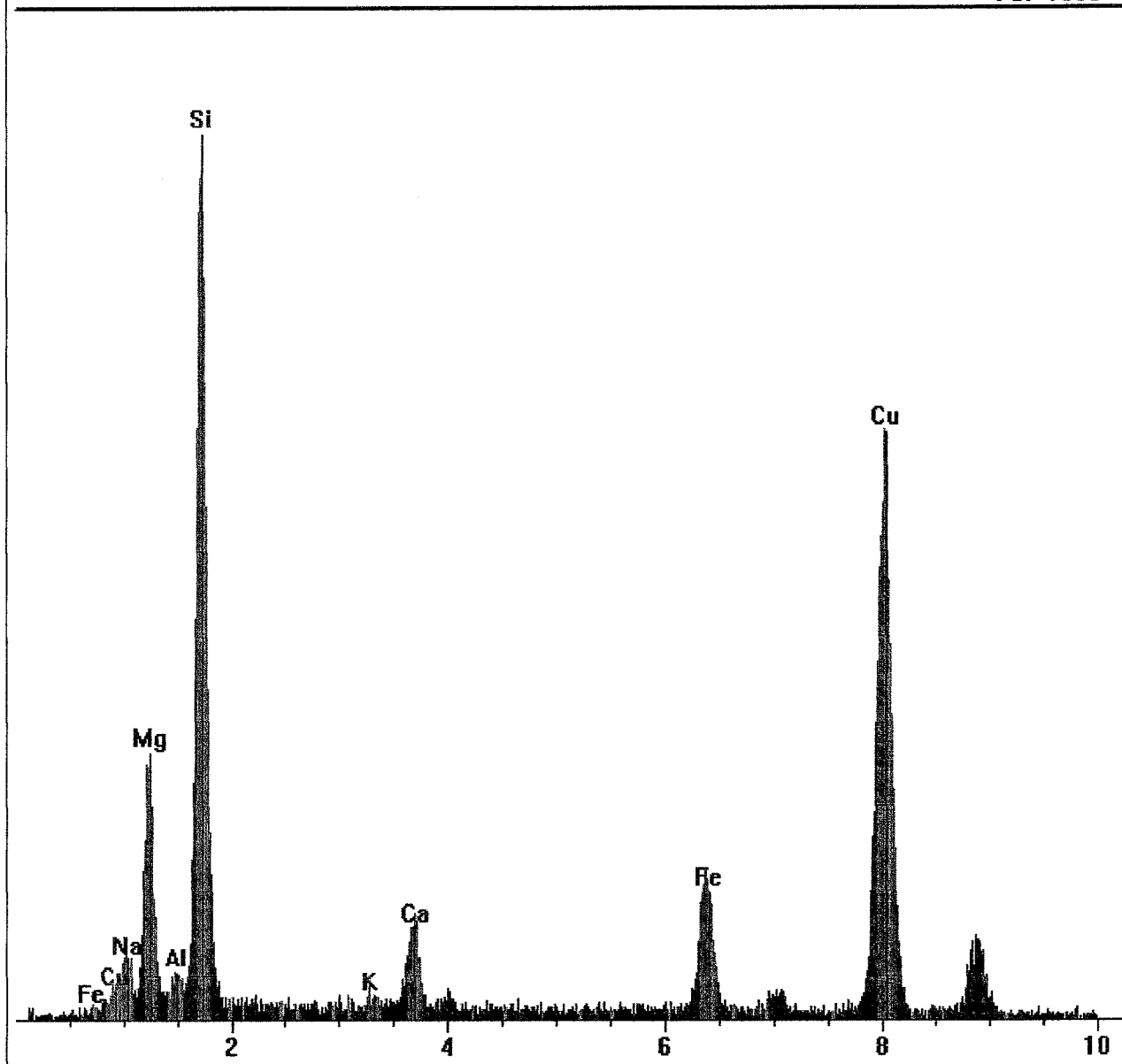

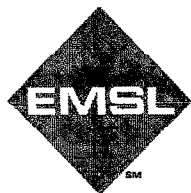

Energy Dispersive X-Ray Analysis Qualitative Spectrum

EMSL ANALYTICAL, INC.

File: F:\Documen...cords\PGT Files\EMSL27-2\EMSL27-2 2013\sw01201\B03.pgt
Collected: September 24, 2013 08:56:43


Live Time: 319.14
Beam Voltage: 20.00


Count Rate: 717
Beam Current: 2.00

Dead Time: 10.38 %
Takeoff Angle: 60.98

■ B03.pgt

FS: 1200

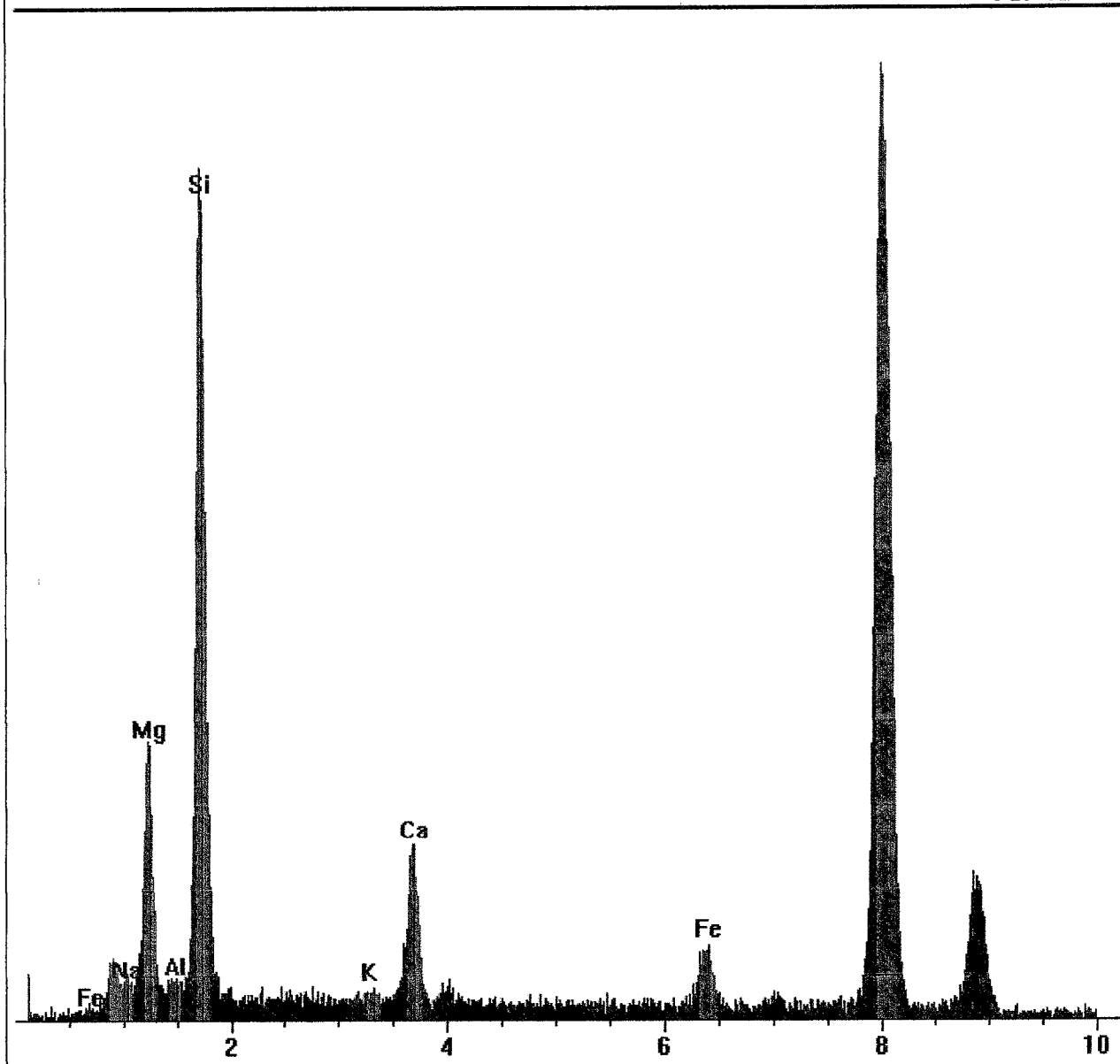

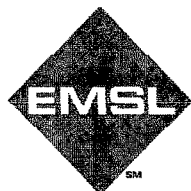

Energy Dispersive X-Ray Analysis

Qualitative Spectrum

EMSL ANALYTICAL, INC.

File: F:\Documen...cords\PGT Files\EMSL27-2\EMSL27-2 2013\sw01201\B04.pgt

Collected: September 24, 2013 09:10:50

Live Time: 100.36

Count Rate: 4453

Dead Time: 40.52 %


Beam Voltage: 20.00

Beam Current: 2.00

Takeoff Angle: 60.98

■ B04.pgt

FS: 2500

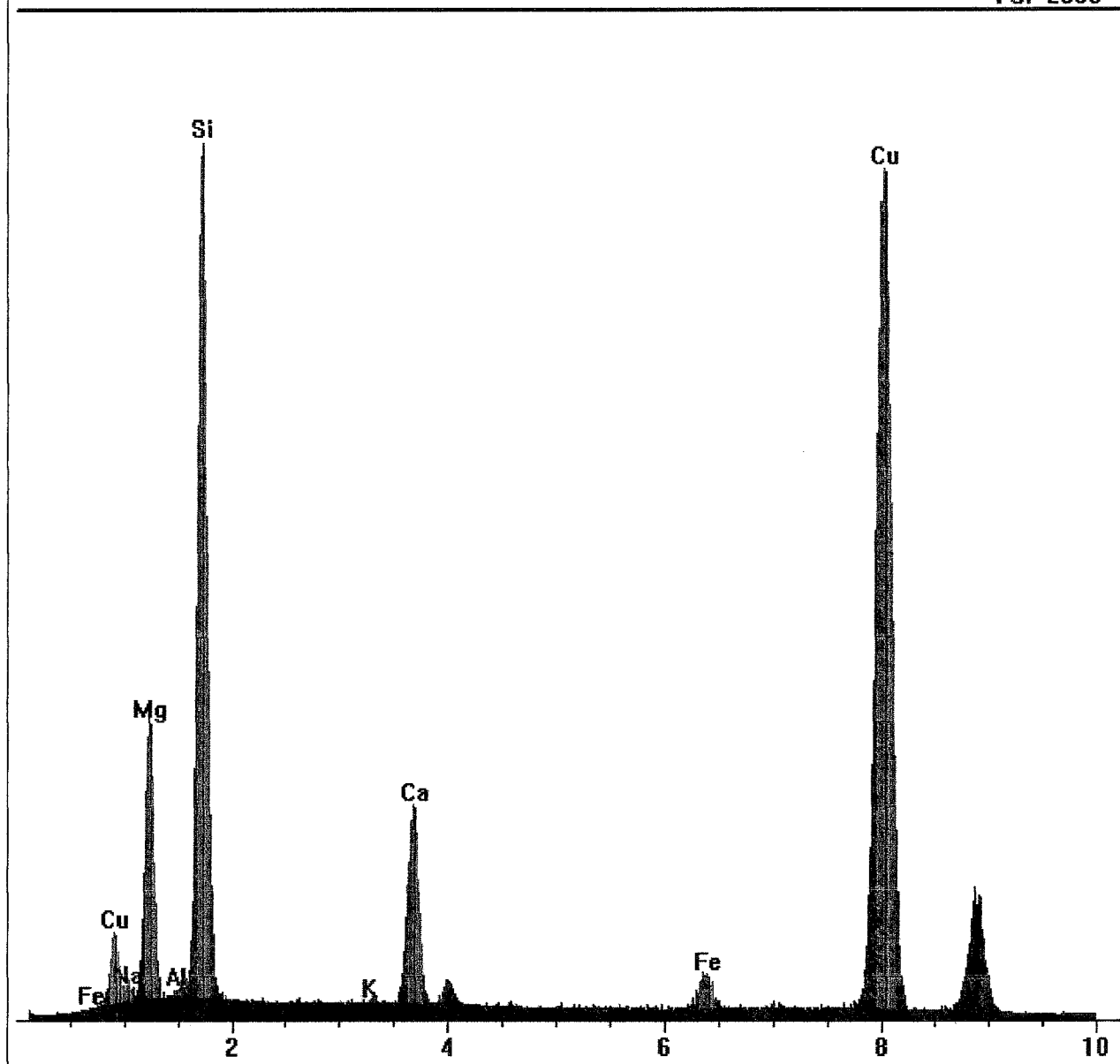



Energy Dispersive X-Ray Analysis

Qualitative Spectrum

EMSL ANALYTICAL, INC.

File: F:\Documen...ords\PGT Files\EMSL27-2\EMSL27-2 2013\sw01201\B05.pgt
Collected: September 24, 2013 09:10:50


Live Time: 50.82
Beam Voltage: 20.00


Count Rate: 3764
Beam Current: 2.00

Dead Time: 35.16 %
Takeoff Angle: 60.98

■ B05.pgt

FS: 1800

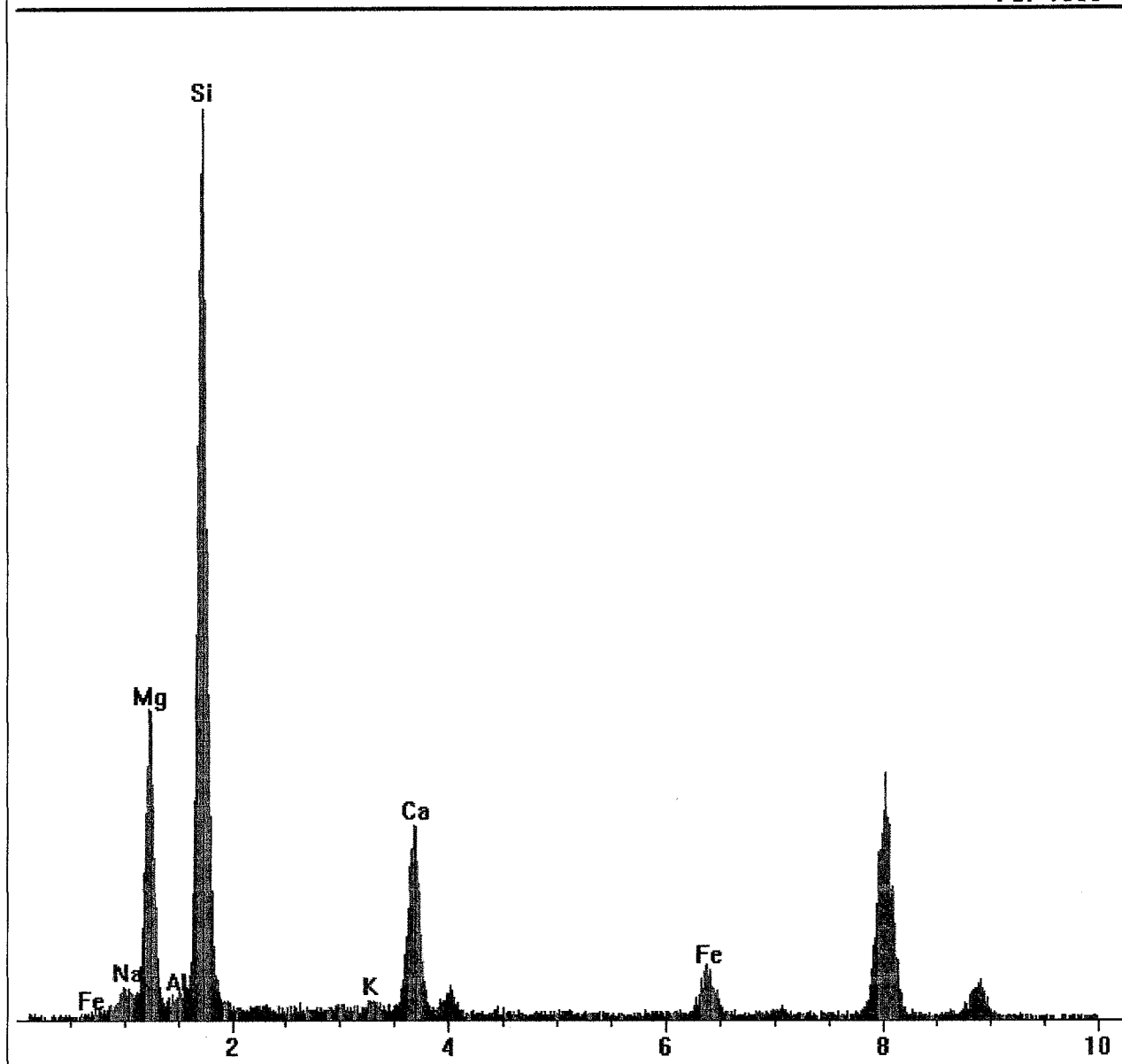

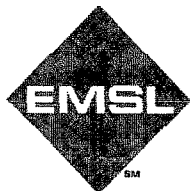

Energy Dispersive X-Ray Analysis Qualitative Spectrum

EMSL ANALYTICAL, INC.

File: F:\Documen...ords\PGT Files\EMSL27-2\EMSL27-2 2013\sw01201\B06.pgt
Collected: September 24, 2013 09:10:50

Live Time: 176.33

Count Rate: 974

Dead Time: 12.41 %


Beam Voltage: 20.00


Beam Current: 2.00

Takeoff Angle: 60.98

■ B06.pgt

FS: 1000

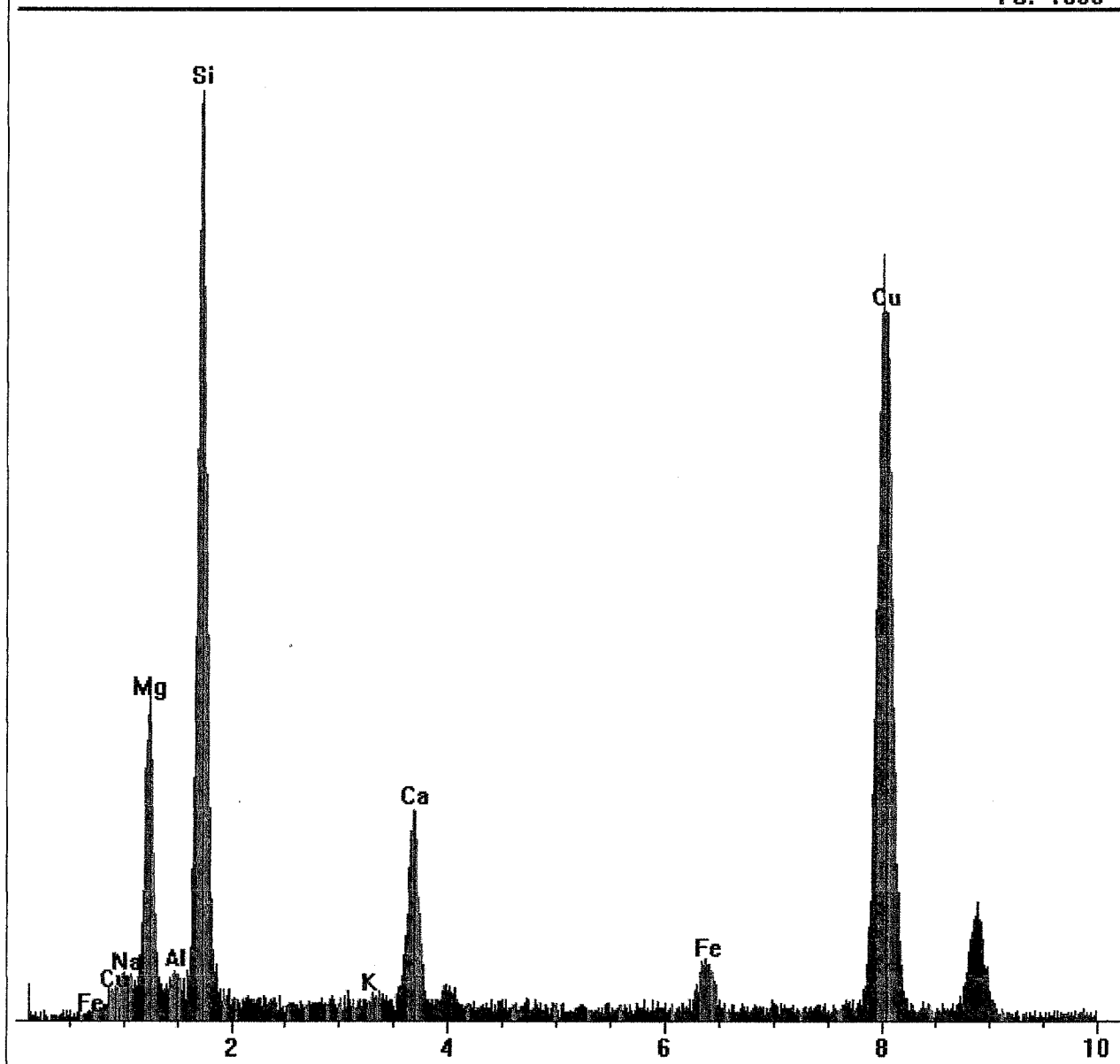

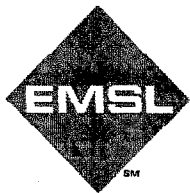

Energy Dispersive X-Ray Analysis Qualitative Spectrum


EMSL ANALYTICAL, INC.


File: F:\Documen...ords\PGT Files\EMSL27-2\EMSL27-2 2013\sw01201\B07.pgt
Collected: September 24, 2013 09:10:50

Live Time: 56.10
Beam Voltage: 20.00

Count Rate: 2114
Beam Current: 2.00

Dead Time: 22.47 %
Takeoff Angle: 60.98

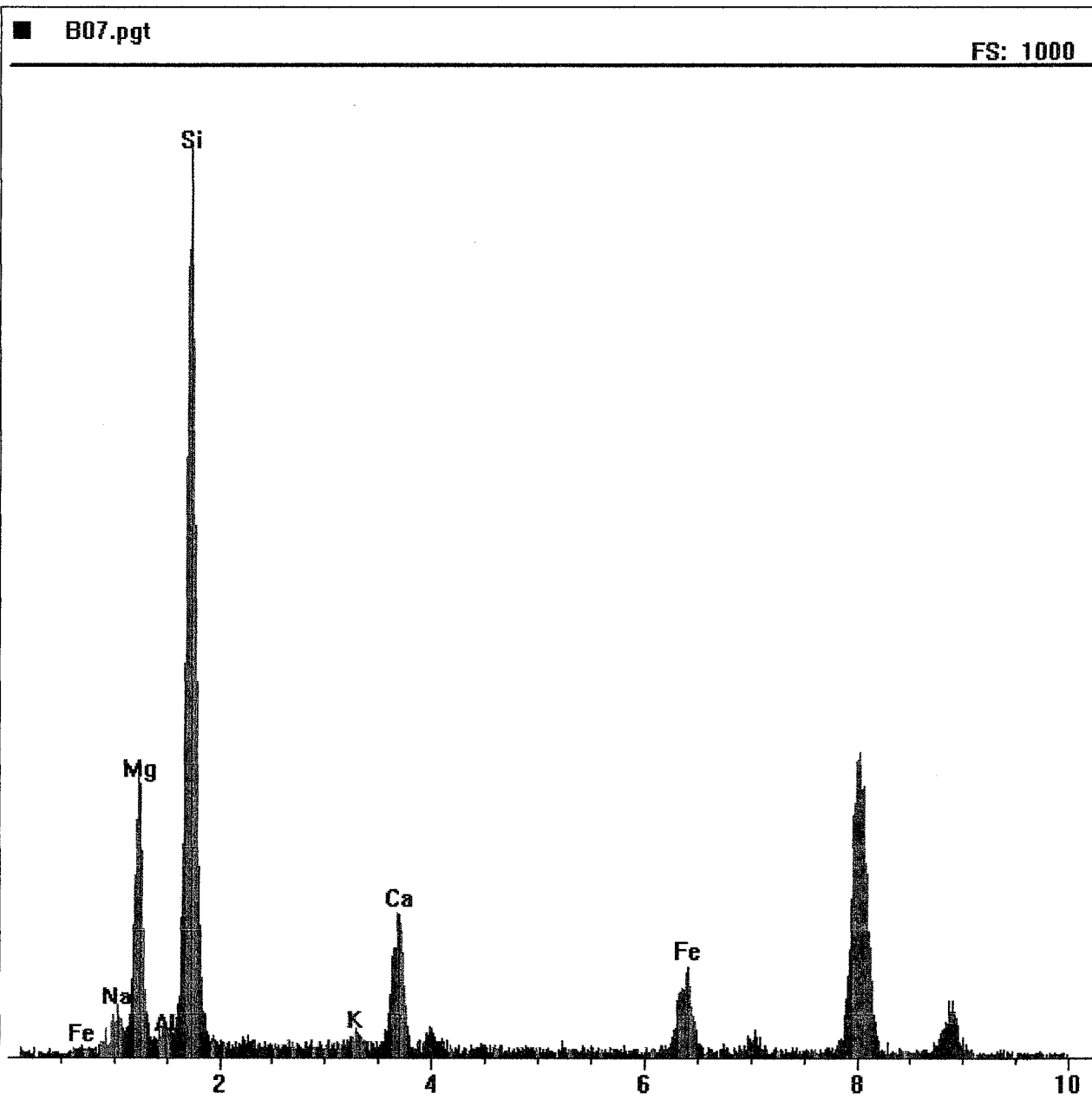

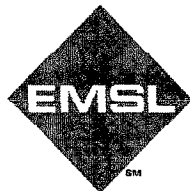

Energy Dispersive X-Ray Analysis

Qualitative Spectrum

EMSL ANALYTICAL, INC.

File: F:\Documen...cords\PGT Files\EMSL27-2\EMSL27-2 2013\sw01201\B08.pgt
Collected: September 24, 2013 09:10:50


Live Time: 56.62
Beam Voltage: 20.00


Count Rate: 3559
Beam Current: 2.00

Dead Time: 33.49 %
Takeoff Angle: 60.98

■ B08.pgt

FS: 1800

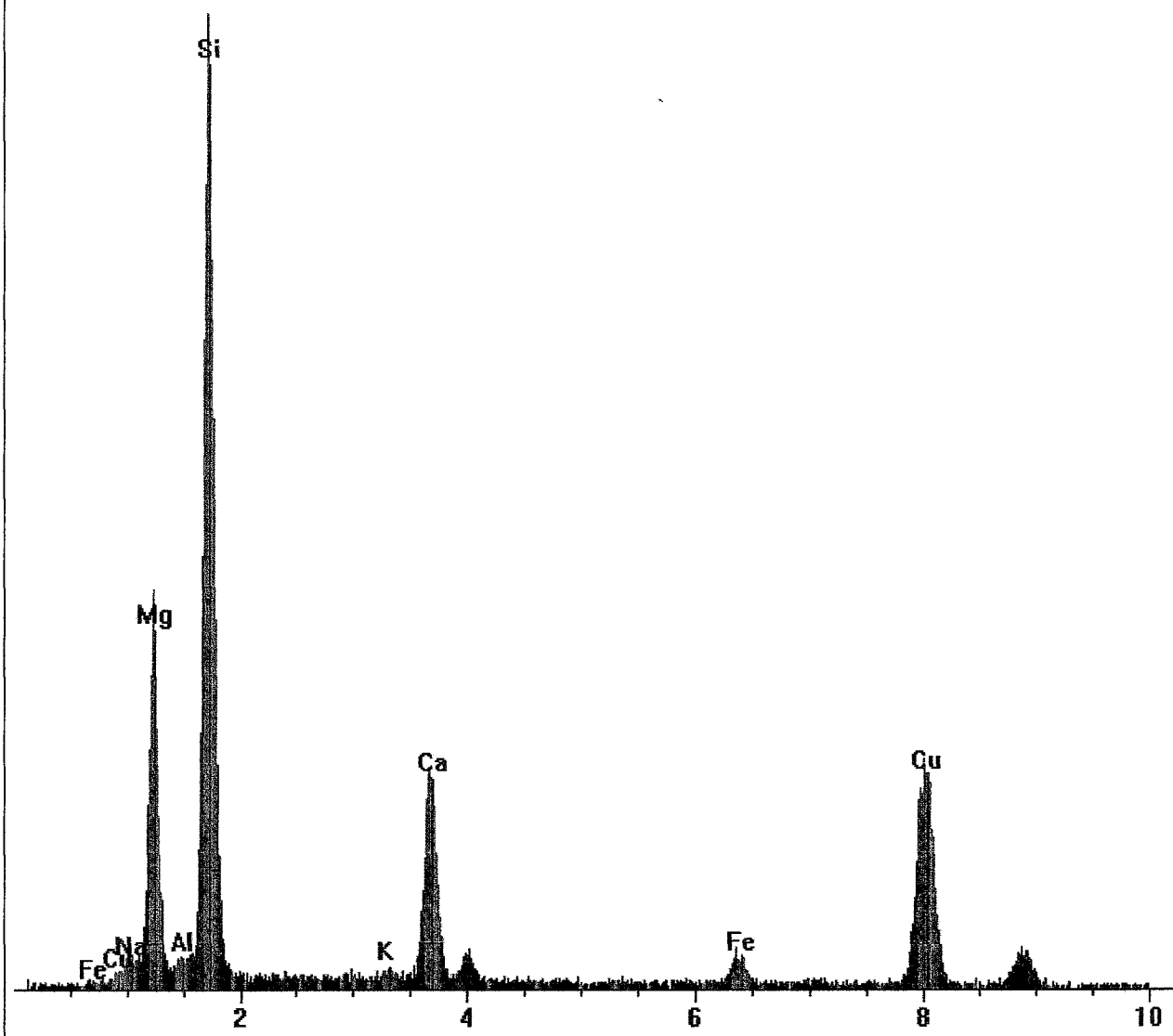

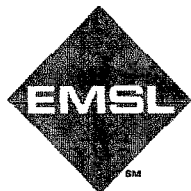

Energy Dispersive X-Ray Analysis

Qualitative Spectrum


EMSL ANALYTICAL, INC.

File: F:\Documen...cords\PGT Files\EMSL27-2\EMSL27-2 2013\sw01201\B09.pgt
Collected: September 24, 2013 09:10:50

Live Time:	393.64	Count Rate:	258	Dead Time:	4.95 %
Beam Voltage:	20.00	Beam Current:	2.00	Takeoff Angle:	60.98

■ B09.pgt

FS: 1000

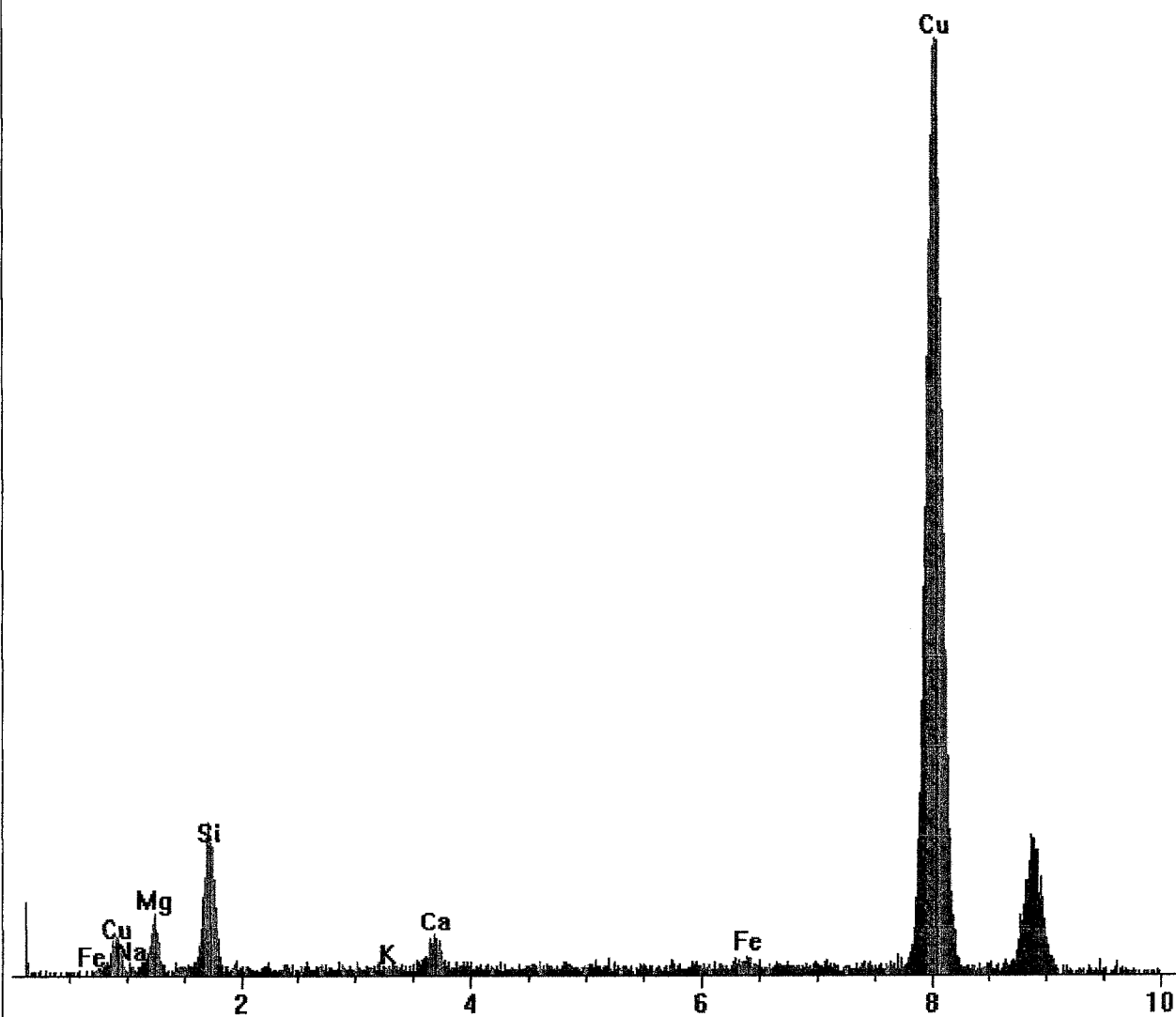



Energy Dispersive X-Ray Analysis

Qualitative Spectrum

EMSL ANALYTICAL, INC.

File: F:\Documen...cords\PGT Files\EMSL27-2\EMSL27-2 2013\sw01201\B10.pgt
Collected: September 24, 2013 09:10:50


Live Time: 121.93
Beam Voltage: 20.00


Count Rate: 1060
Beam Current: 2.00

Dead Time: 14.01 %
Takeoff Angle: 60.98

■ B10.pgt

FS: 1000

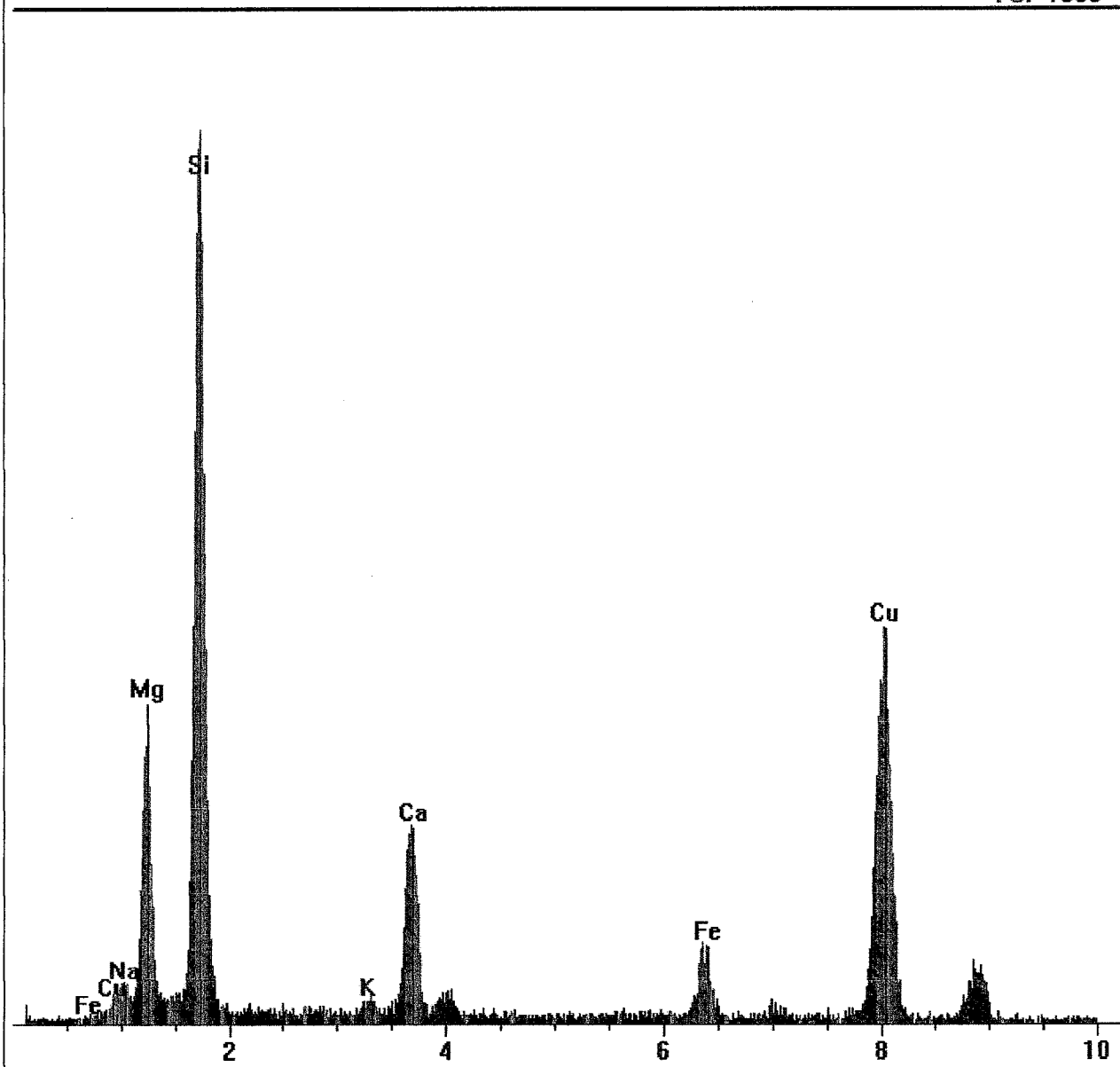

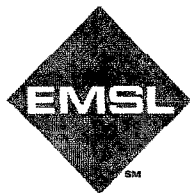

Energy Dispersive X-Ray Analysis Qualitative Spectrum


EMSL ANALYTICAL, INC.


File: F:\Documen...cords\PGT Files\EMSL27-2\EMSL27-2 2013\SW01231\D01.pgt
Collected: September 24, 2013 09:40:13

Live Time: 162.75
Beam Voltage: 20.00

Count Rate: 3645
Beam Current: 2.00

Dead Time: 34.74 %
Takeoff Angle: 60.98

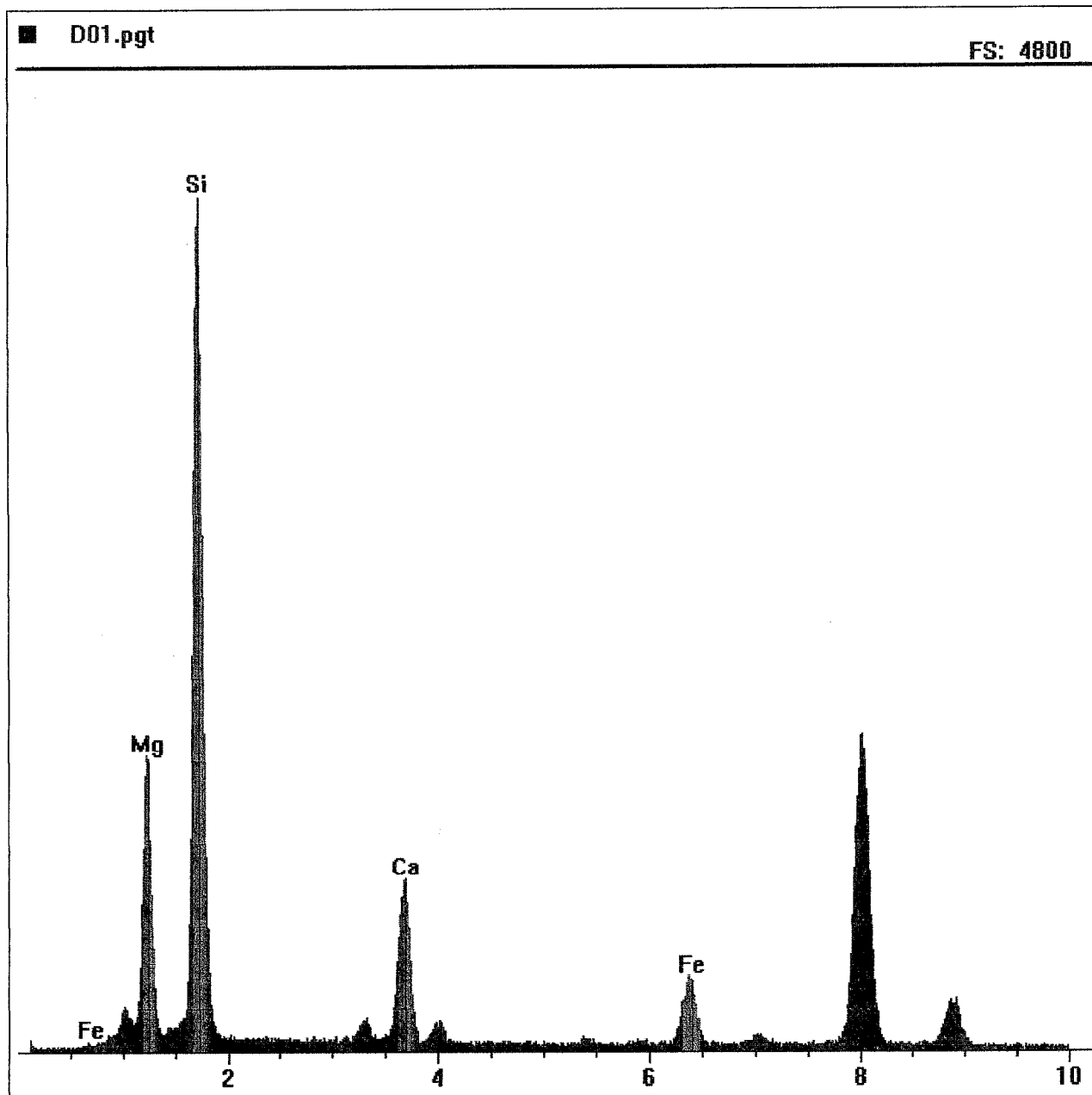

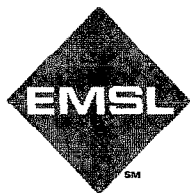

Energy Dispersive X-Ray Analysis

Qualitative Spectrum

EMSL ANALYTICAL, INC.

File: F:\Documen...cords\PGT Files\EMSL27-2\EMSL27-2 2013\SW01231\D02.pgt

Collected: September 24, 2013 09:40:13

Live Time: 281.14

Count Rate: 2478

Dead Time: 26.69 %


Beam Voltage: 20.00


Beam Current: 2.00

Takeoff Angle: 60.98

■ D02.pgt

FS: 4800

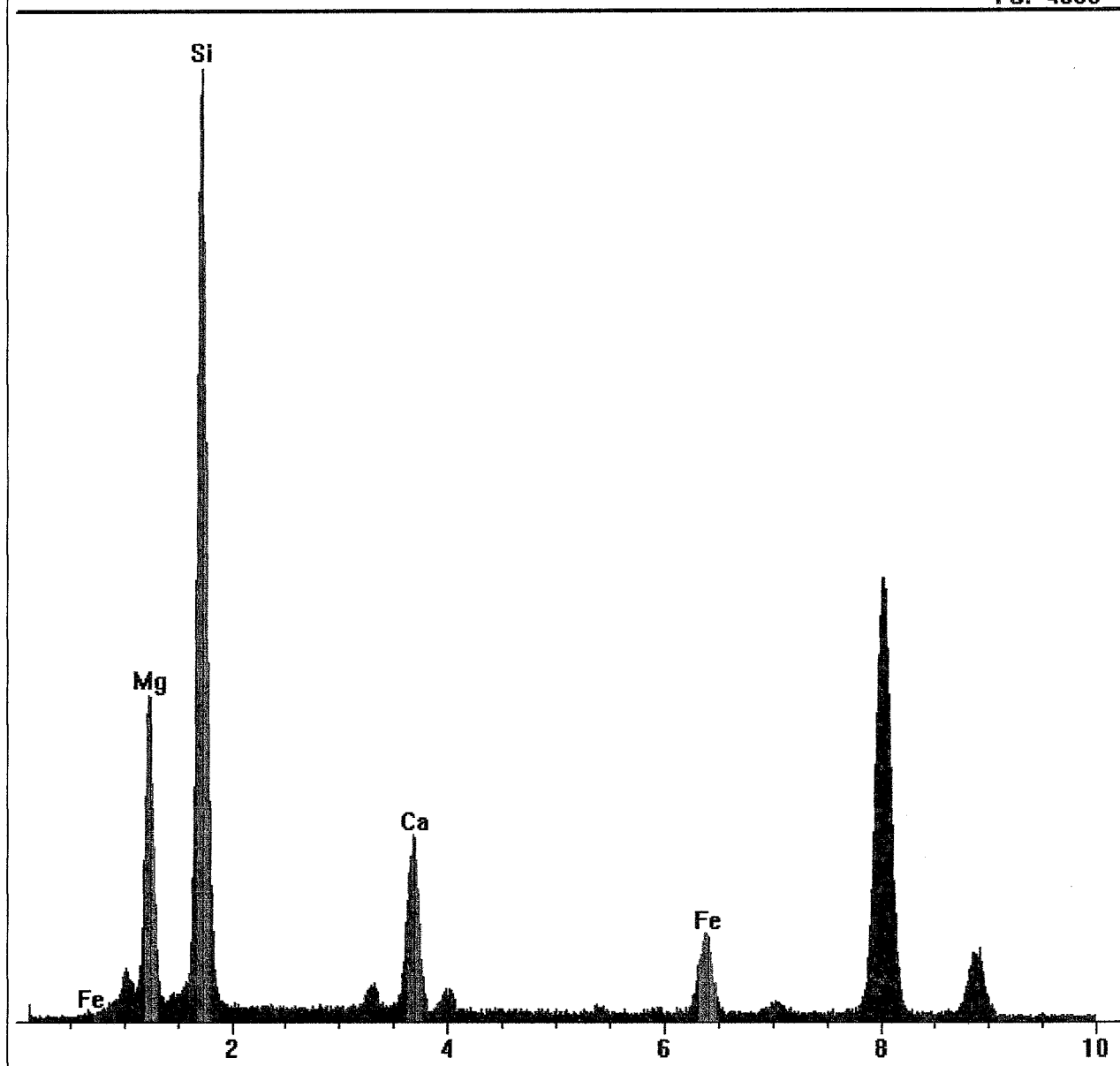

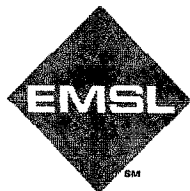


Energy Dispersive X-Ray Analysis


Qualitative Spectrum

EMSL ANALYTICAL, INC.

File: F:\Documen...ords\PGT Files\EMSL27-2\EMSL27-2 2013\SW01231\D03.pgt
Collected: September 24, 2013 09:48:02

Live Time:	272.30	Count Rate:	424	Dead Time:	6.98 %
Beam Voltage:	20.00	Beam Current:	2.00	Takeoff Angle:	60.98

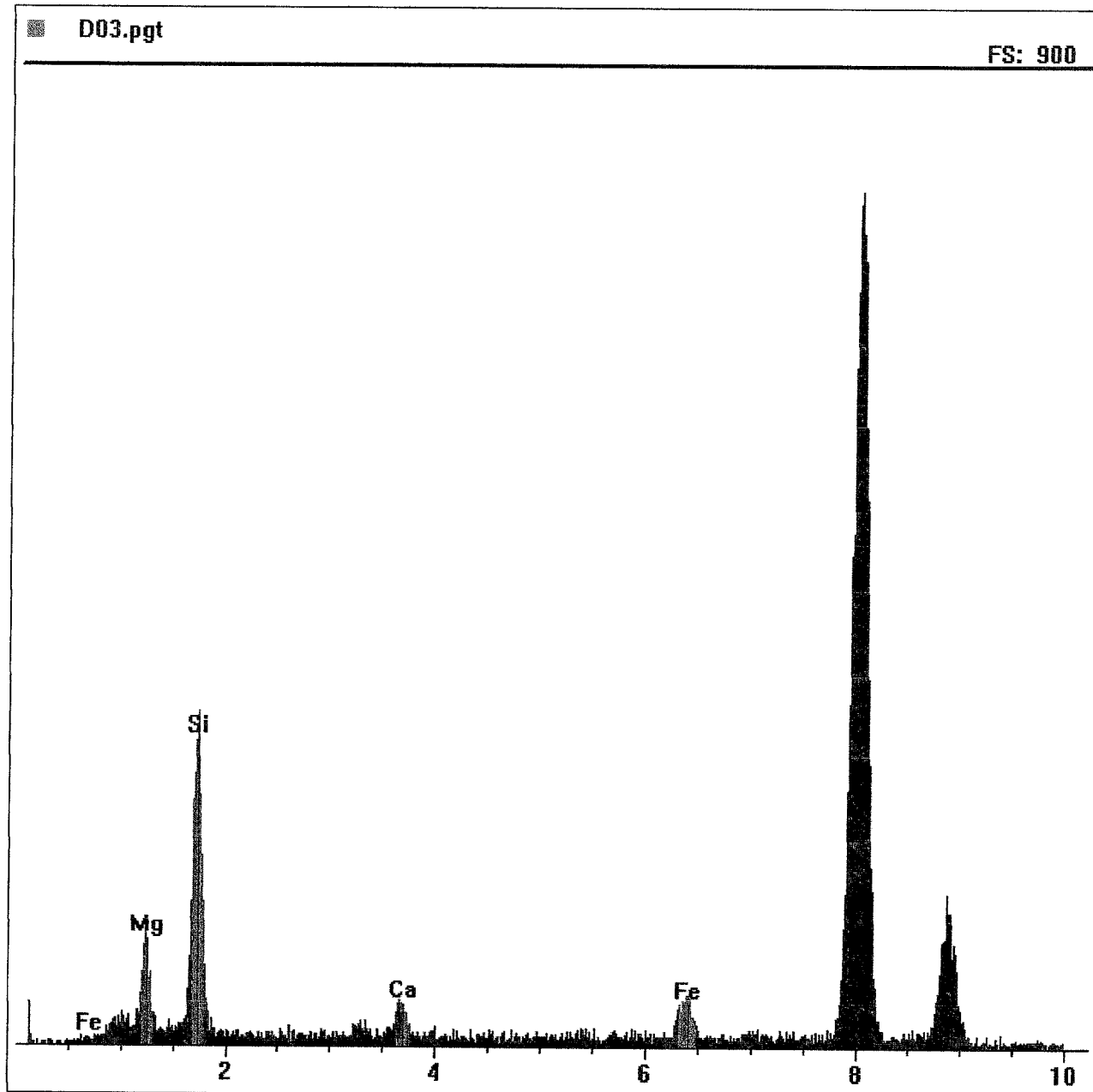

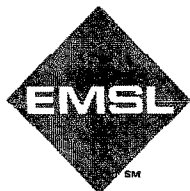

Energy Dispersive X-Ray Analysis Qualitative Spectrum

EMSL ANALYTICAL, INC.

File: F:\Documen...cords\PGT Files\EMSL27-2\EMSL27-2 2013\SW01231\D04.pgt

Collected: September 24, 2013 09:48:02

Live Time: 150.91

Count Rate: 2011

Dead Time: 23.05 %


Beam Voltage: 20.00


Beam Current: 2.00

Takeoff Angle: 60.98

■ D04.pgt

FS: 2000

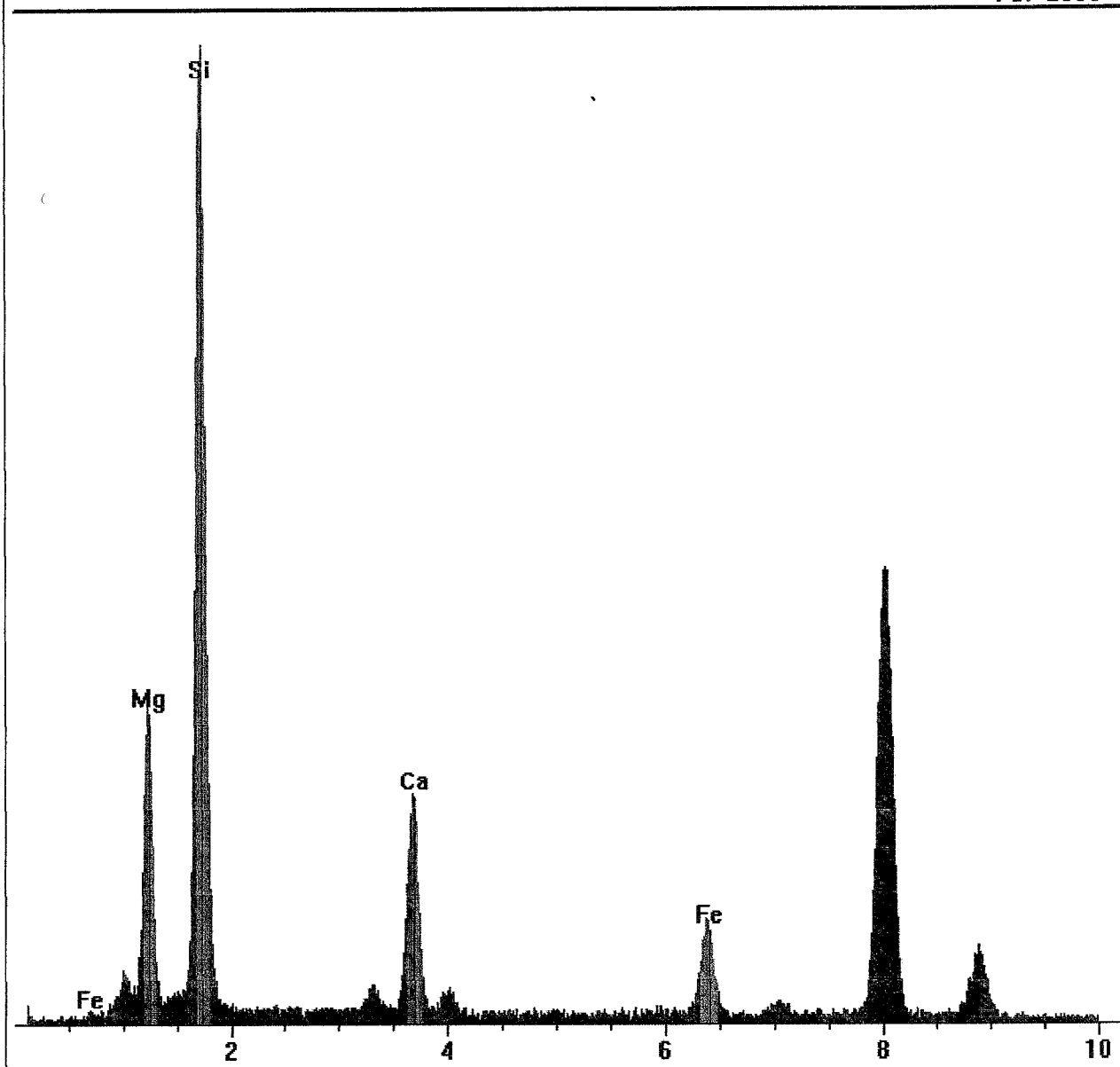

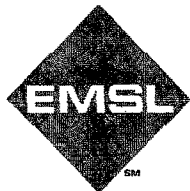

Energy Dispersive X-Ray Analysis Qualitative Spectrum


EMSL ANALYTICAL, INC.


File: F:\Documen...ords\PGT Files\EMSL27-2\EMSL27-2 2013\SW01231\D05.pgt
Collected: September 24, 2013 09:48:02

Live Time: 117.84
Beam Voltage: 20.00

Count Rate: 1863
Beam Current: 2.00

Dead Time: 22.72 %
Takeoff Angle: 60.98

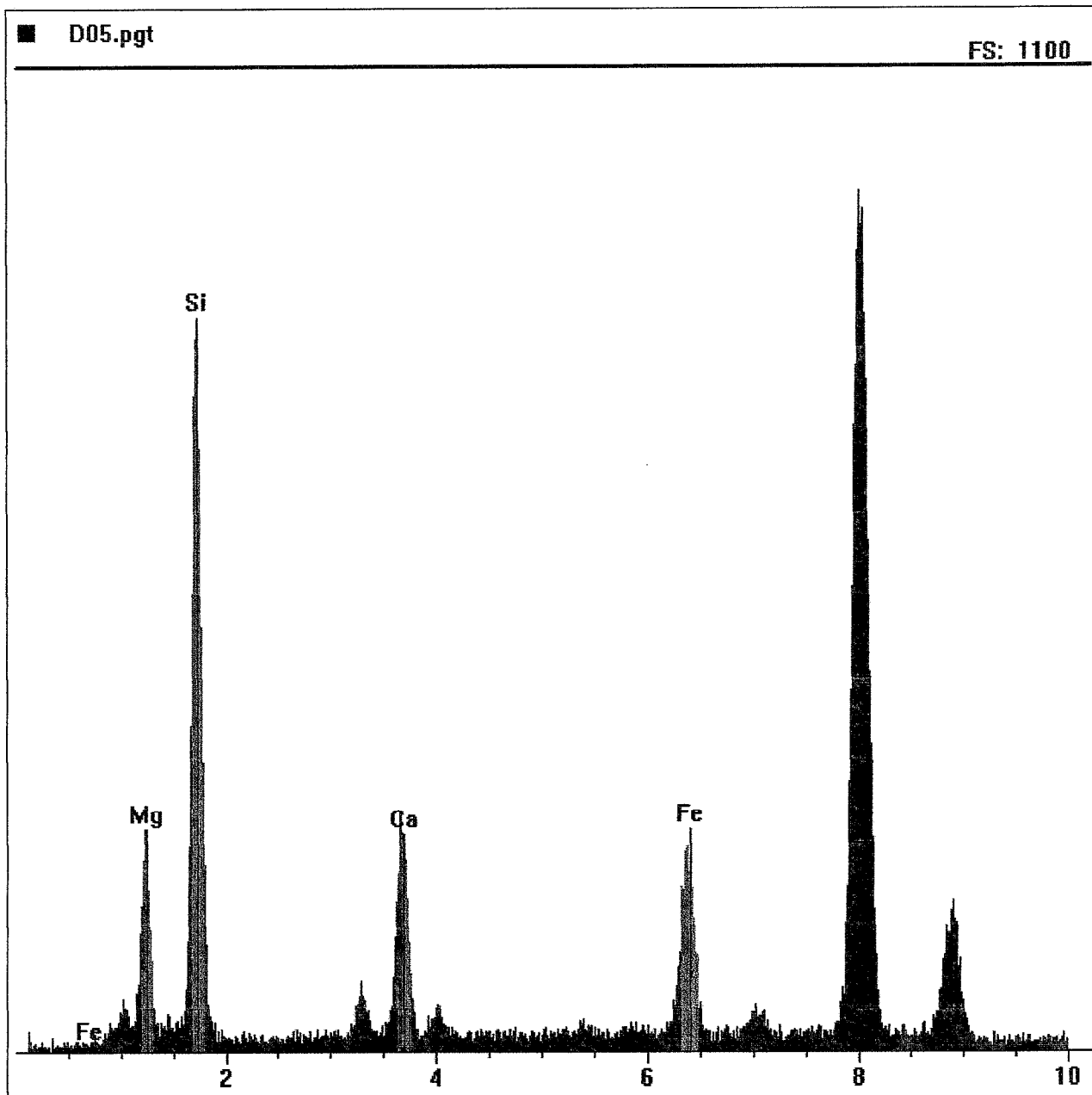

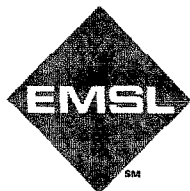

Energy Dispersive X-Ray Analysis

Qualitative Spectrum

EMSL ANALYTICAL, INC.

File: F:\Documen...ords\PGT Files\EMSL27-2\EMSL27-2 2013\SW01231\D06.pgt
Collected: September 24, 2013 09:48:02


Live Time: 74.78
Beam Voltage: 20.00


Count Rate: 1883
Beam Current: 2.00

Dead Time: 20.75 %
Takeoff Angle: 60.98

■ D06.pgt

FS: 1100

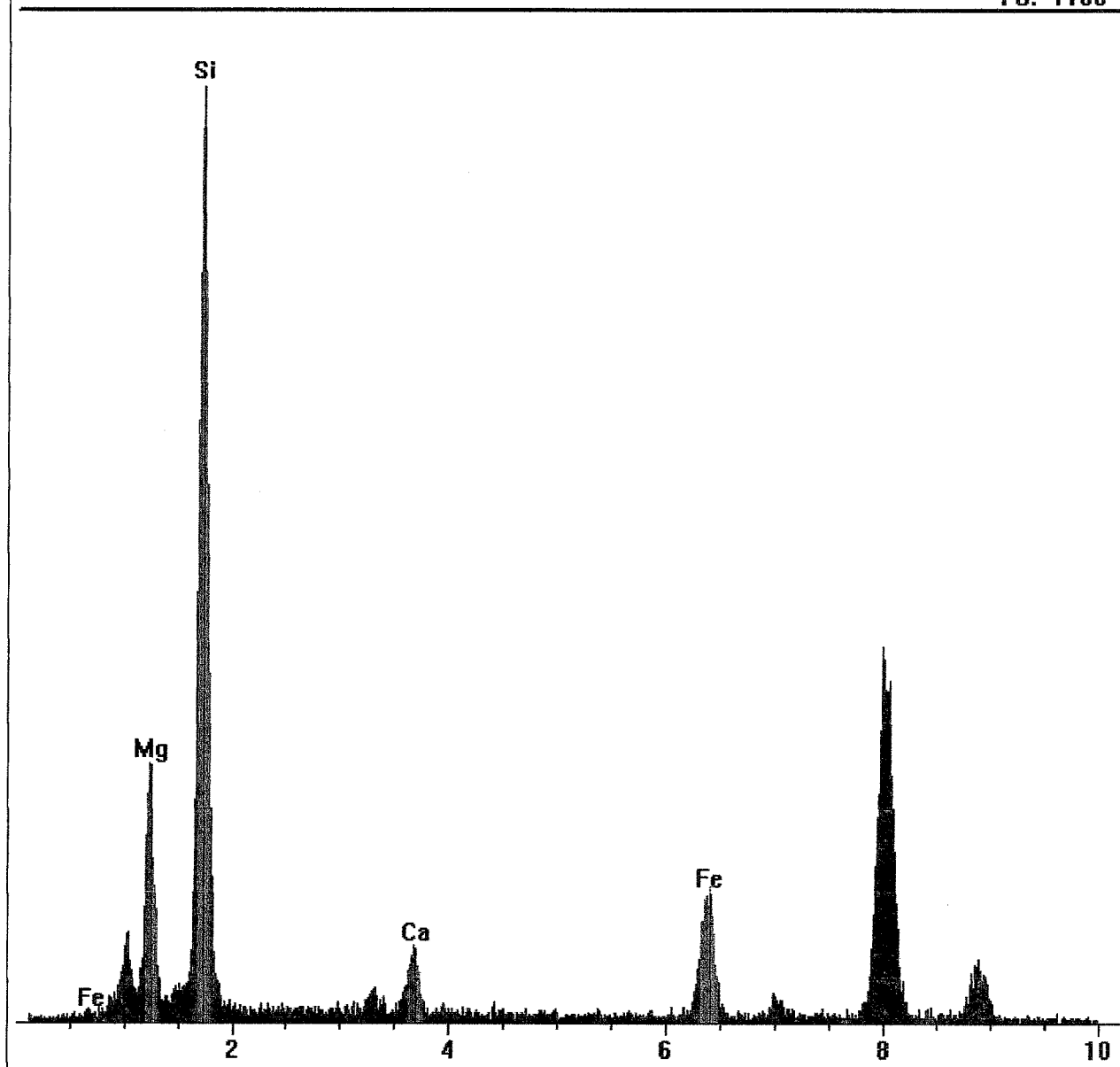

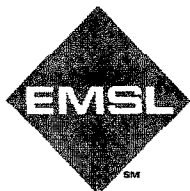

Energy Dispersive X-Ray Analysis Qualitative Spectrum

EMSL ANALYTICAL, INC.

File: F:\Documen...cords\PGT Files\EMSL27-2\EMSL27-2 2013\SW01231\D07.pgt
Collected: September 24, 2013 09:48:02


Live Time: 119.62
Beam Voltage: 20.00


Count Rate: 1284
Beam Current: 2.00

Dead Time: 15.29 %
Takeoff Angle: 60.98

■ D07.pgt

FS: 1000

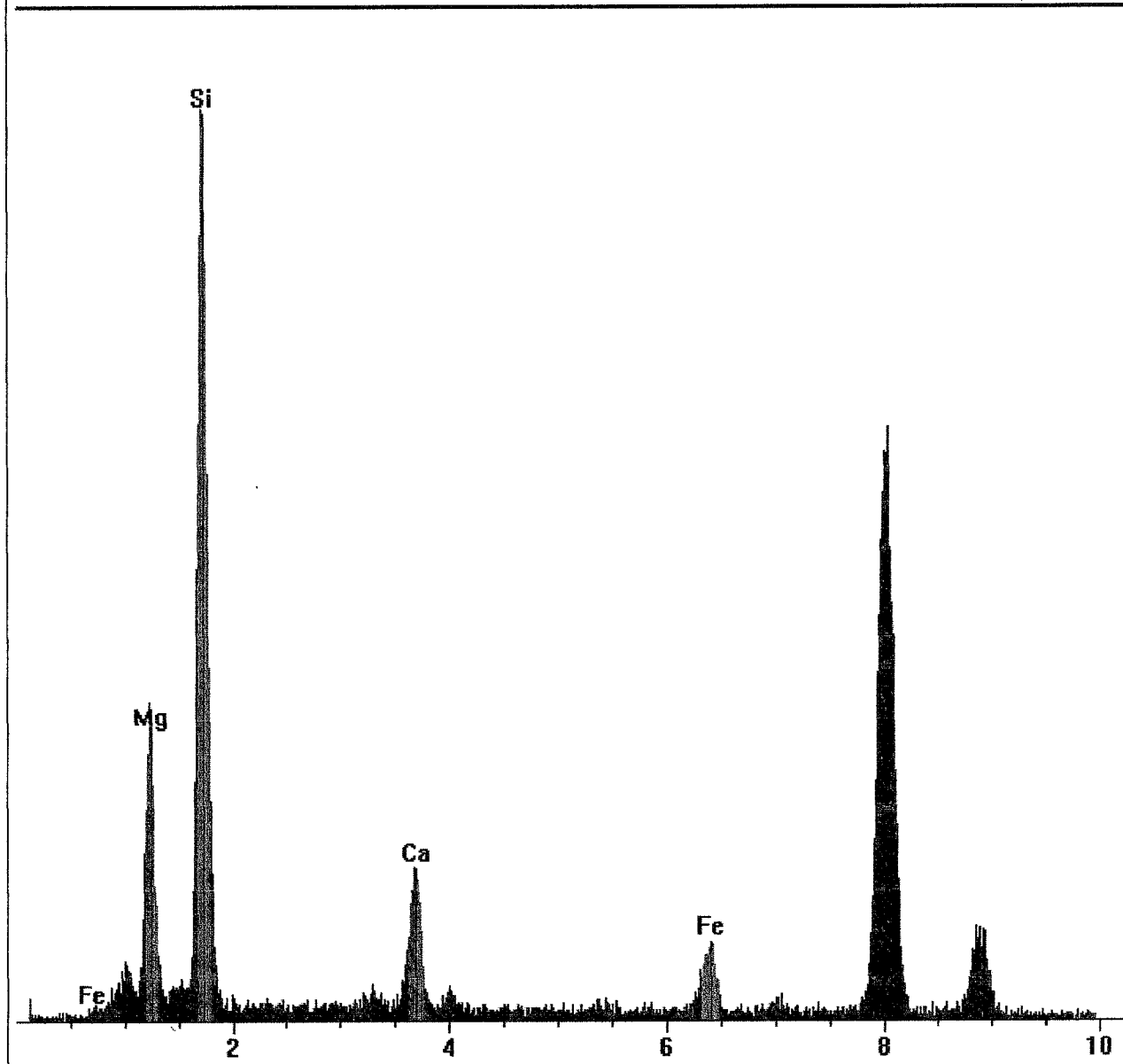

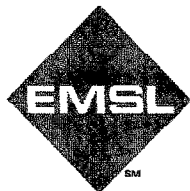

Energy Dispersive X-Ray Analysis Qualitative Spectrum


EMSL ANALYTICAL, INC.

File: F:\Documen...ords\PGT Files\EMSL27-2\EMSL27-2 2013\SW01231\D08.pgt
Collected: September 24, 2013 09:48:02

Live Time: 173.64
Beam Voltage: 20.00

Count Rate: 1409
Beam Current: 2.00

Dead Time: 16.40 %
Takeoff Angle: 60.98

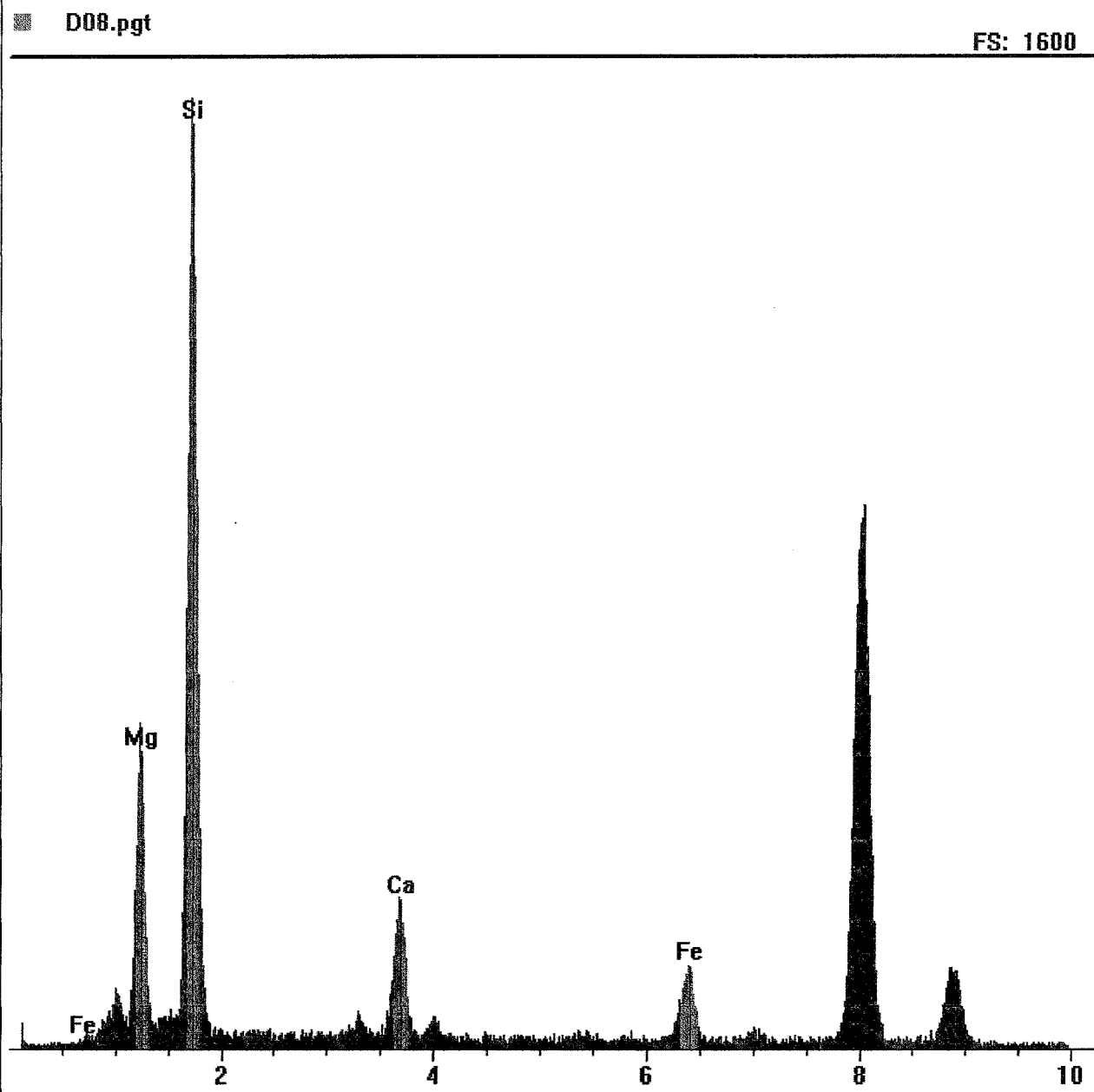



Energy Dispersive X-Ray Analysis

Qualitative Spectrum

EMSL ANALYTICAL, INC.

File: F:\Documen...ords\PGT Files\EMSL27-2\EMSL27-2 2013\SW01231\D09.pgt
Collected: September 24, 2013 09:48:02

Live Time: 173.64

Count Rate: 1409

Dead Time: 16.40 %


Beam Voltage: 20.00


Beam Current: 2.00

Takeoff Angle: 60.98

■ D09.pgt

FS: 1600

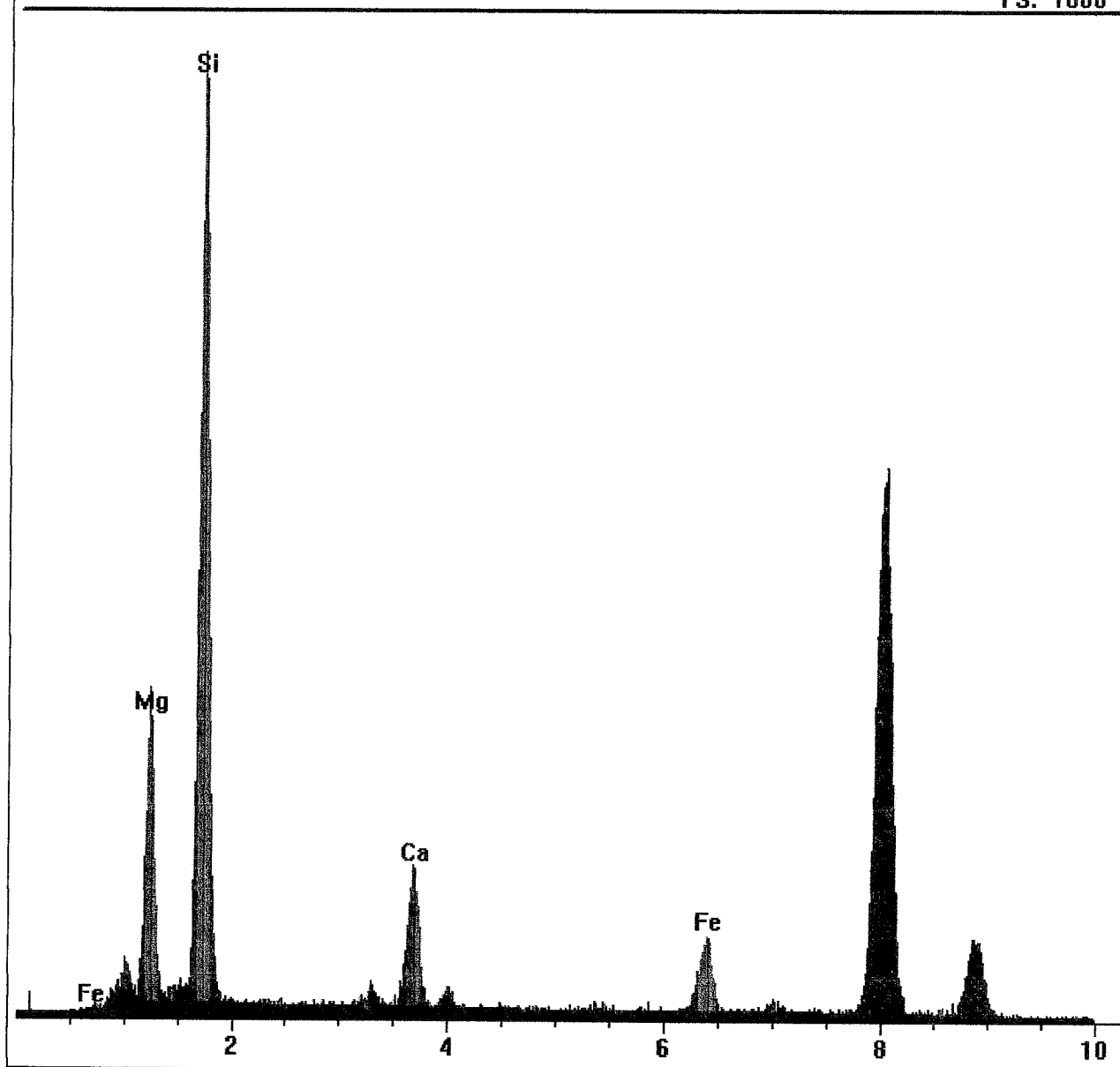

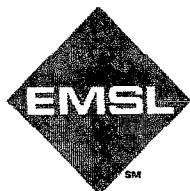

Energy Dispersive X-Ray Analysis

Qualitative Spectrum

EMSL ANALYTICAL, INC.

File: F:\Documen...ords\PGT Files\EMSL27-2\EMSL27-2 2013\SW01231\D10.pgt

Collected: September 24, 2013 09:48:02

Live Time: 209.32


Count Rate: 931


Dead Time: 12.20 %

Beam Voltage: 20.00

Beam Current: 2.00

Takeoff Angle: 60.98

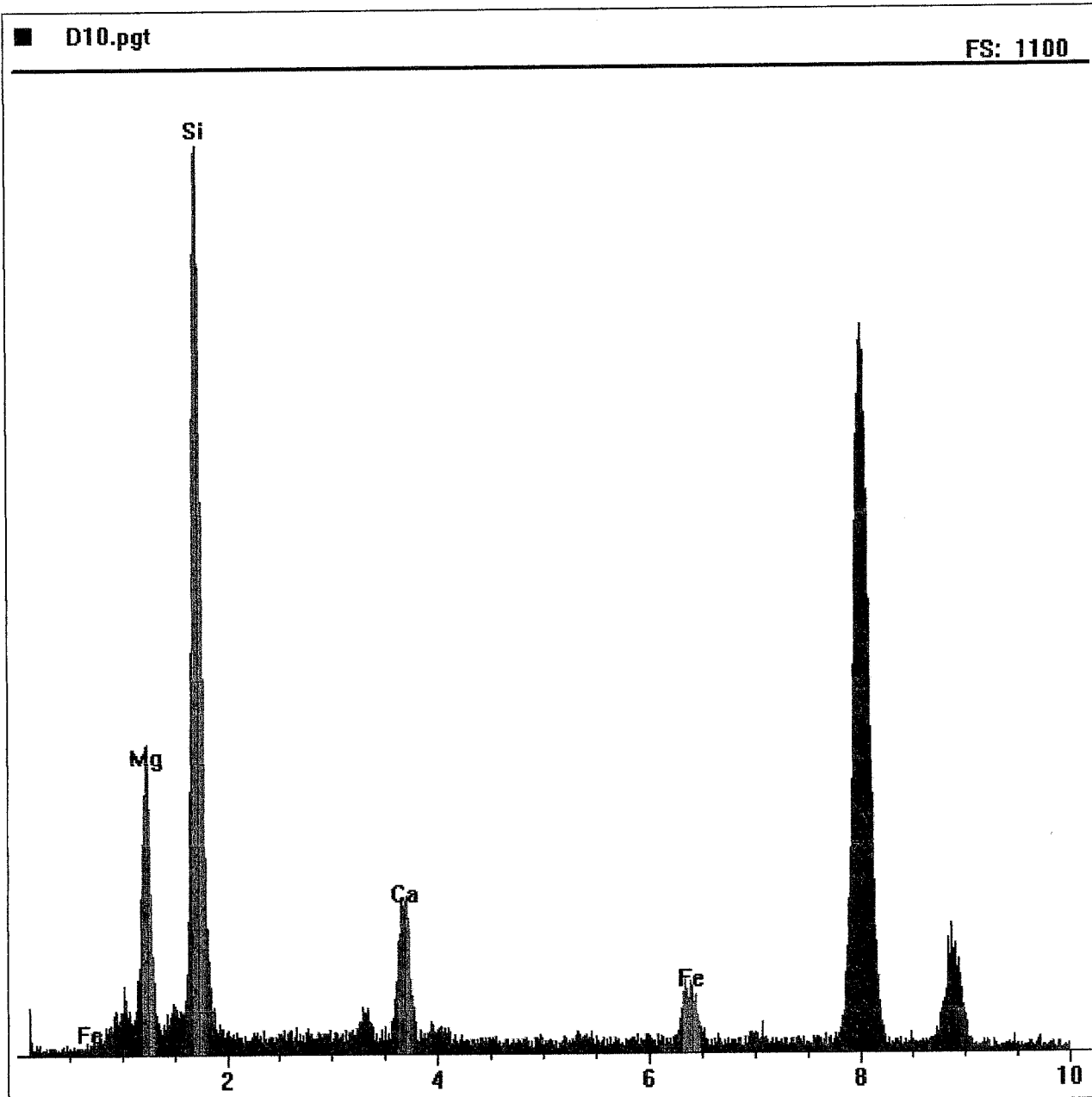

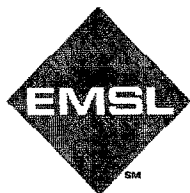

Energy Dispersive X-Ray Analysis

Qualitative Spectrum


EMSL ANALYTICAL, INC.


File: F:\Documen...ecords\PGT Files\EMSL27-2\EMSL27-2 2013\SW21281\E1.pgt
Collected: September 23, 2013 11:08:12

Live Time: 484.60
Beam Voltage: 20.00

Count Rate: 284
Beam Current: 2.00

Dead Time: 5.14 %
Takeoff Angle: 60.98

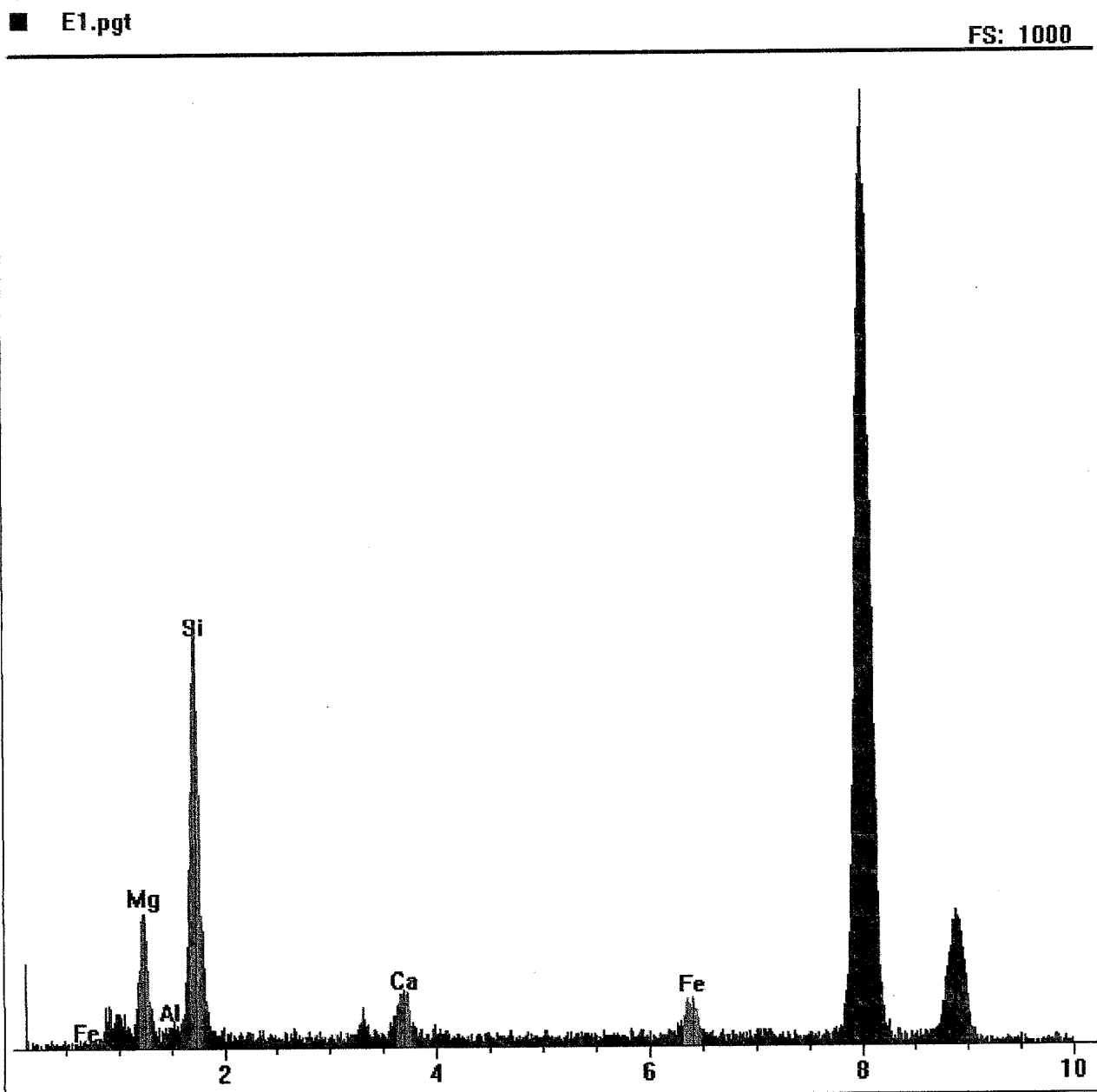

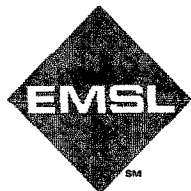

Energy Dispersive X-Ray Analysis

Qualitative Spectrum


EMSL ANALYTICAL, INC.


File: F:\Documen...ecords\PGT Files\EMSL27-2\EMSL27-2 2013\SW21281\E2.pgt
Collected: September 23, 2013 11:08:12

Live Time: 910.26
Beam Voltage: 20.00

Count Rate: 177
Beam Current: 2.00

Dead Time: 3.93 %
Takeoff Angle: 60.98

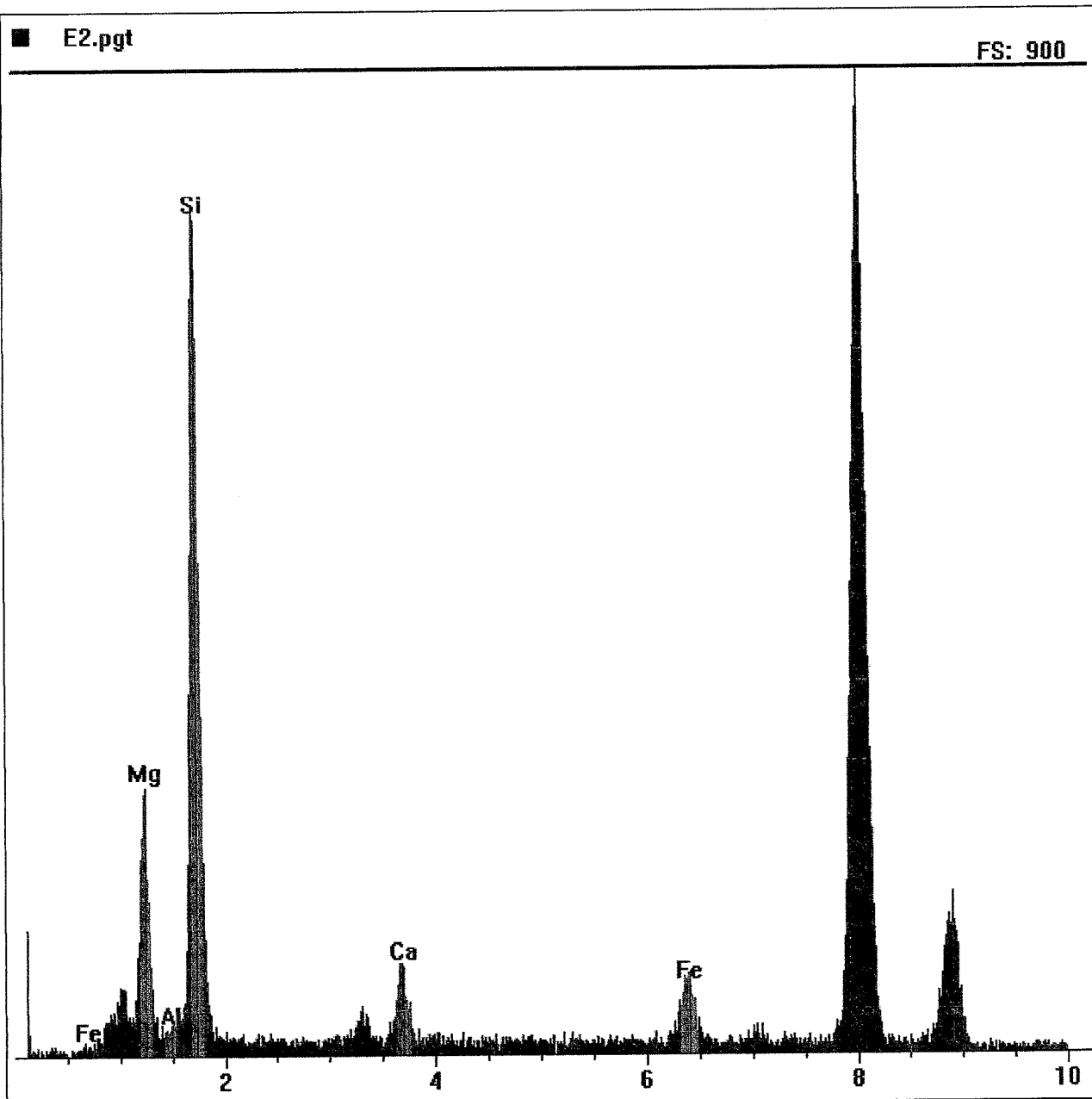

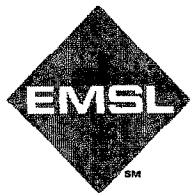

Energy Dispersive X-Ray Analysis Qualitative Spectrum

EMSL ANALYTICAL, INC.

File: F:\Documen...ecords\PGT Files\EMSL27-2\EMSL27-2 2013\SW21281\E3.pgt
Collected: September 23, 2013 11:08:12


Live Time: 583.64
Beam Voltage: 20.00


Count Rate: 607
Beam Current: 2.00

Dead Time: 8.91 %
Takeoff Angle: 60.98

■ E3.pgt

FS: 2000

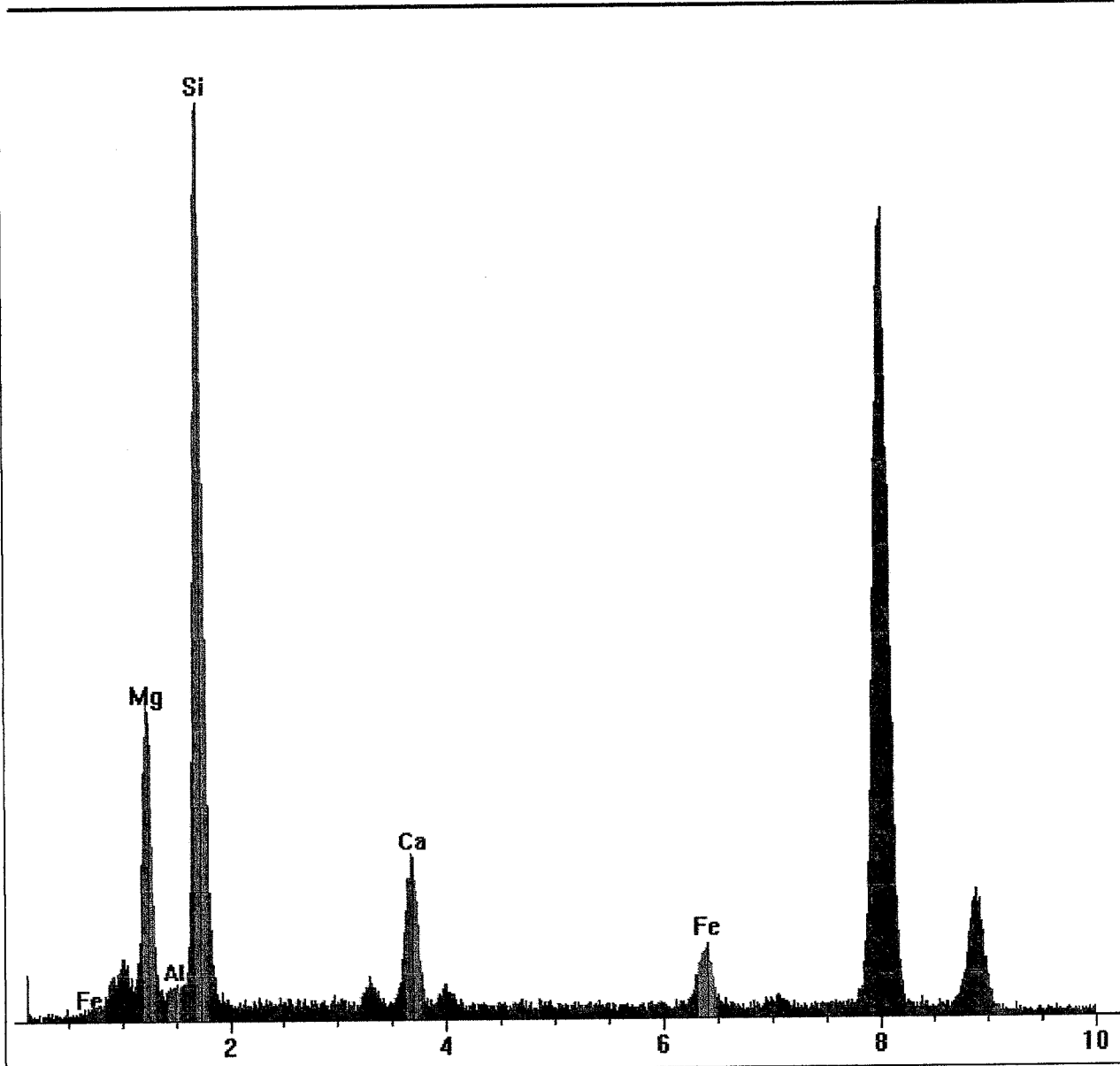

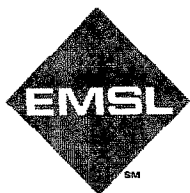

Energy Dispersive X-Ray Analysis

Qualitative Spectrum


EMSL ANALYTICAL, INC.


File: F:\Documen...ecords\PGT Files\EMSL27-2\EMSL27-2 2013\SW21281\E4.pgt
Collected: September 24, 2013 09:48:02

Live Time: 63.22
Beam Voltage: 20.00

Count Rate: 3066
Beam Current: 2.00

Dead Time: 30.24 %
Takeoff Angle: 60.98

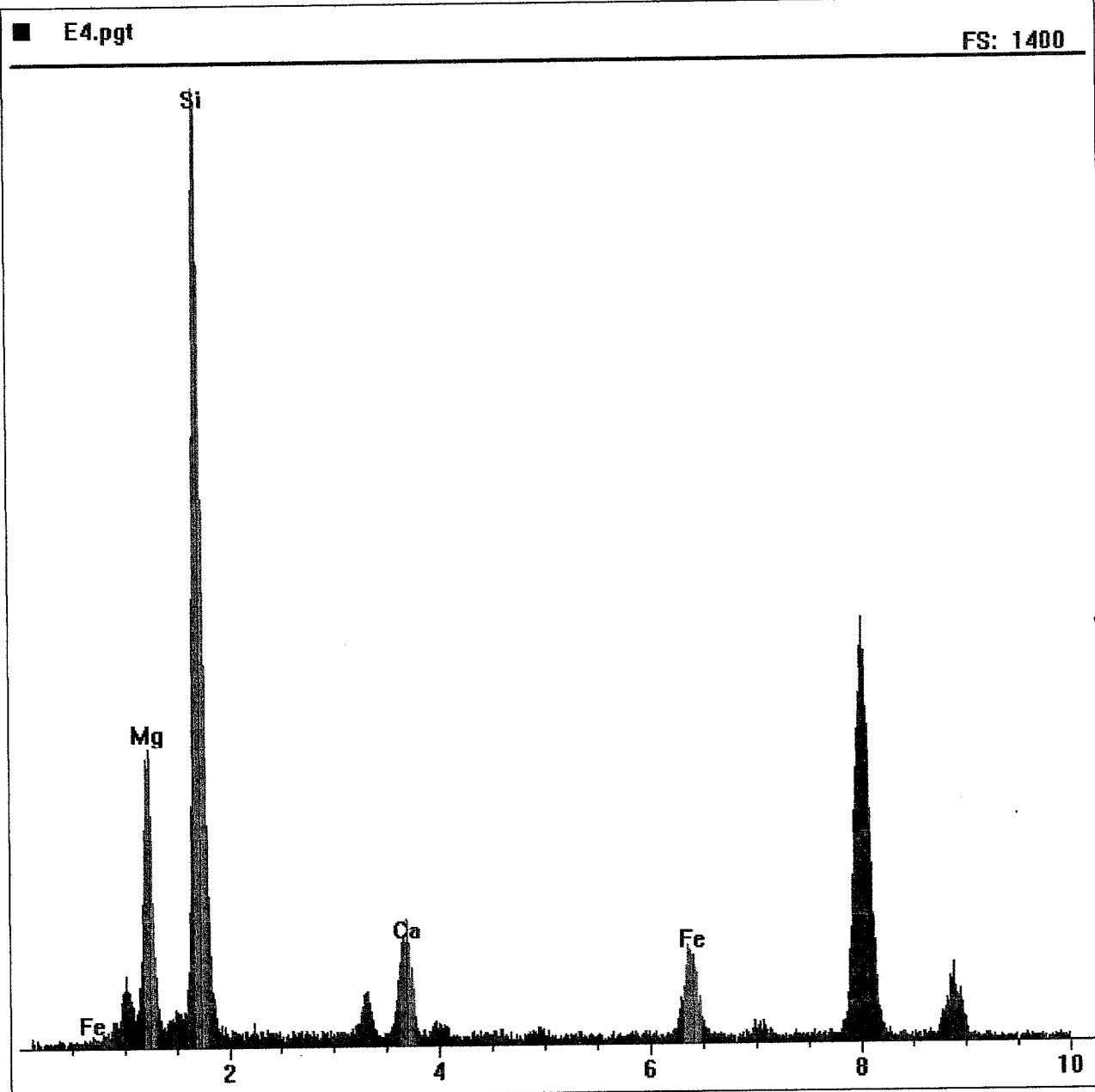

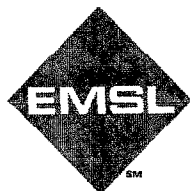

Energy Dispersive X-Ray Analysis

Qualitative Spectrum


EMSL ANALYTICAL, INC.

File: F:\Documen...ecords\PGT Files\EMSL27-2\EMSL27-2 2013\SW21281\E5.pgt
Collected: September 24, 2013 09:48:02

Live Time: 230.37
Beam Voltage: 20.00

Count Rate: 1644
Beam Current: 2.00

Dead Time: 19.08 %
Takeoff Angle: 60.98

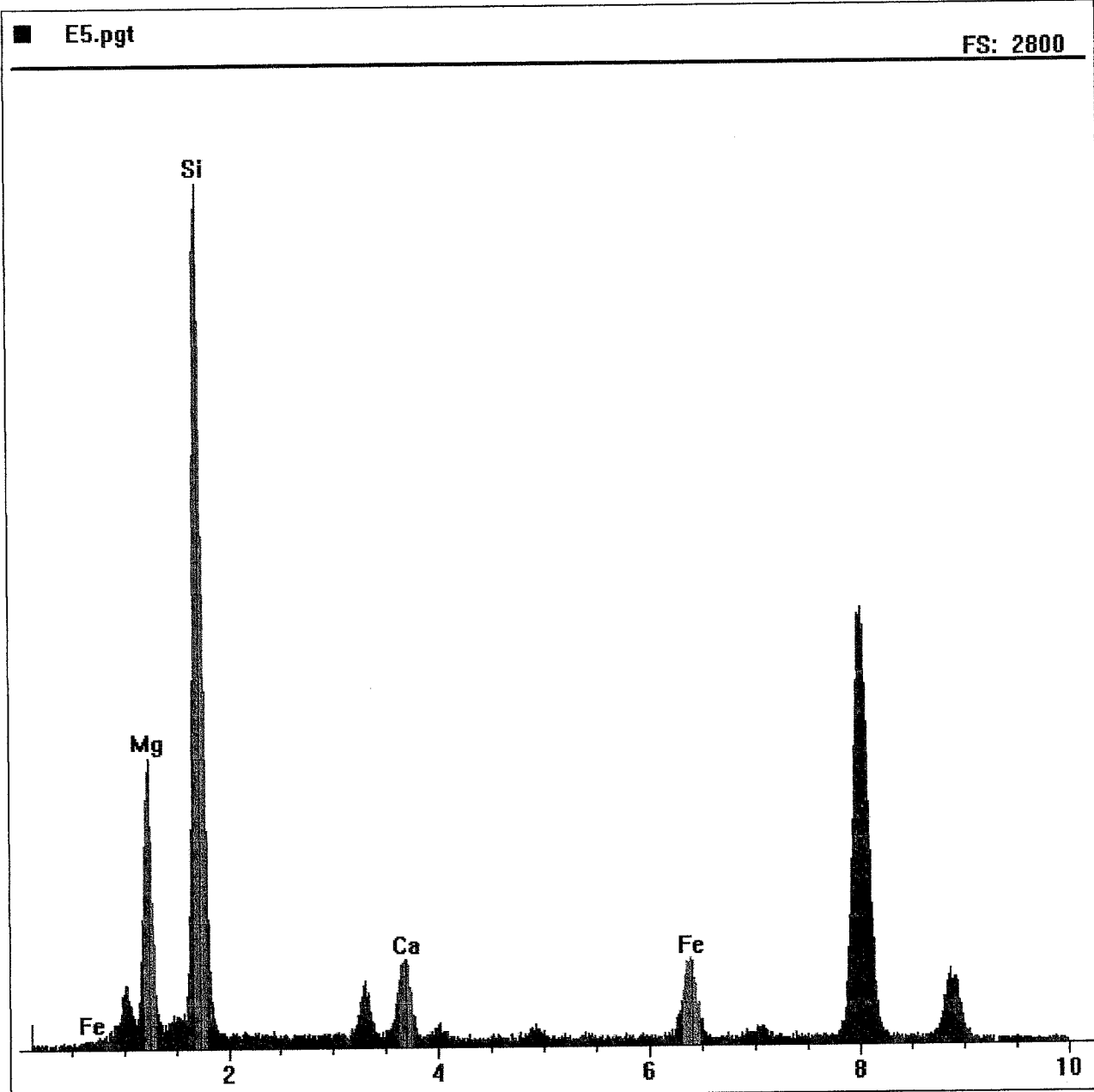



Energy Dispersive X-Ray Analysis

Qualitative Spectrum

EMSL ANALYTICAL, INC.

File: F:\Documen...ecords\PGT Files\EMSL27-2\EMSL27-2 2013\SW21281\E6.pgt
Collected: September 23, 2013 11:08:12


Live Time: 484.60
Beam Voltage: 20.00


Count Rate: 284
Beam Current: 2.00

Dead Time: 5.14 %
Takeoff Angle: 60.98

■ E6.pgt

FS: 1000

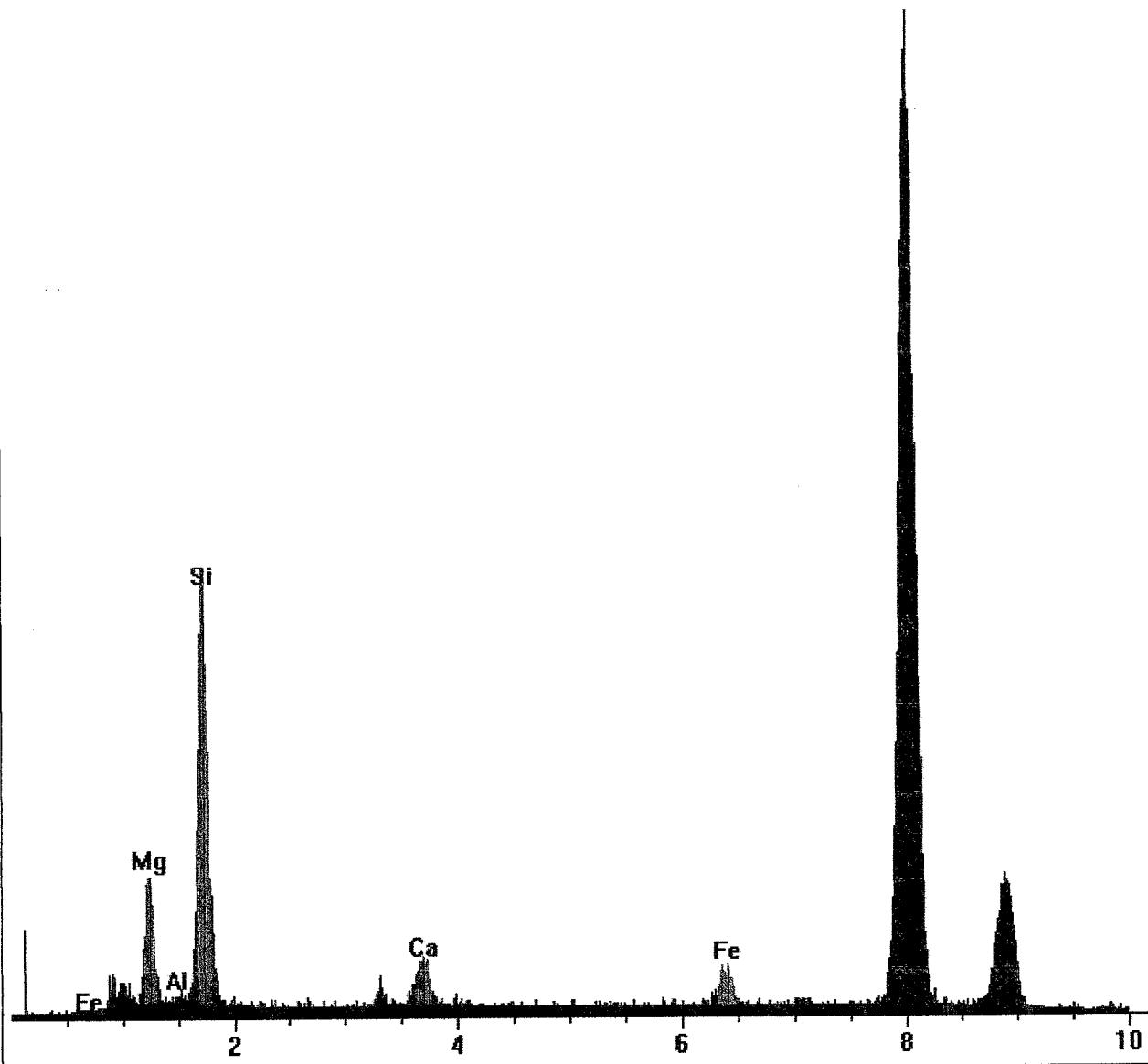

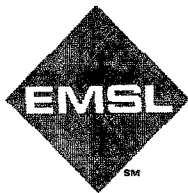

Energy Dispersive X-Ray Analysis Qualitative Spectrum


EMSL ANALYTICAL, INC.


File: F:\Documen...ecords\PGT Files\EMSL27-2\EMSL27-2 2013\SW21281\E7.pgt
Collected: September 24, 2013 09:48:02

Live Time: 59.20
Beam Voltage: 20.00

Count Rate: 2353
Beam Current: 2.00

Dead Time: 24.49 %
Takeoff Angle: 60.98

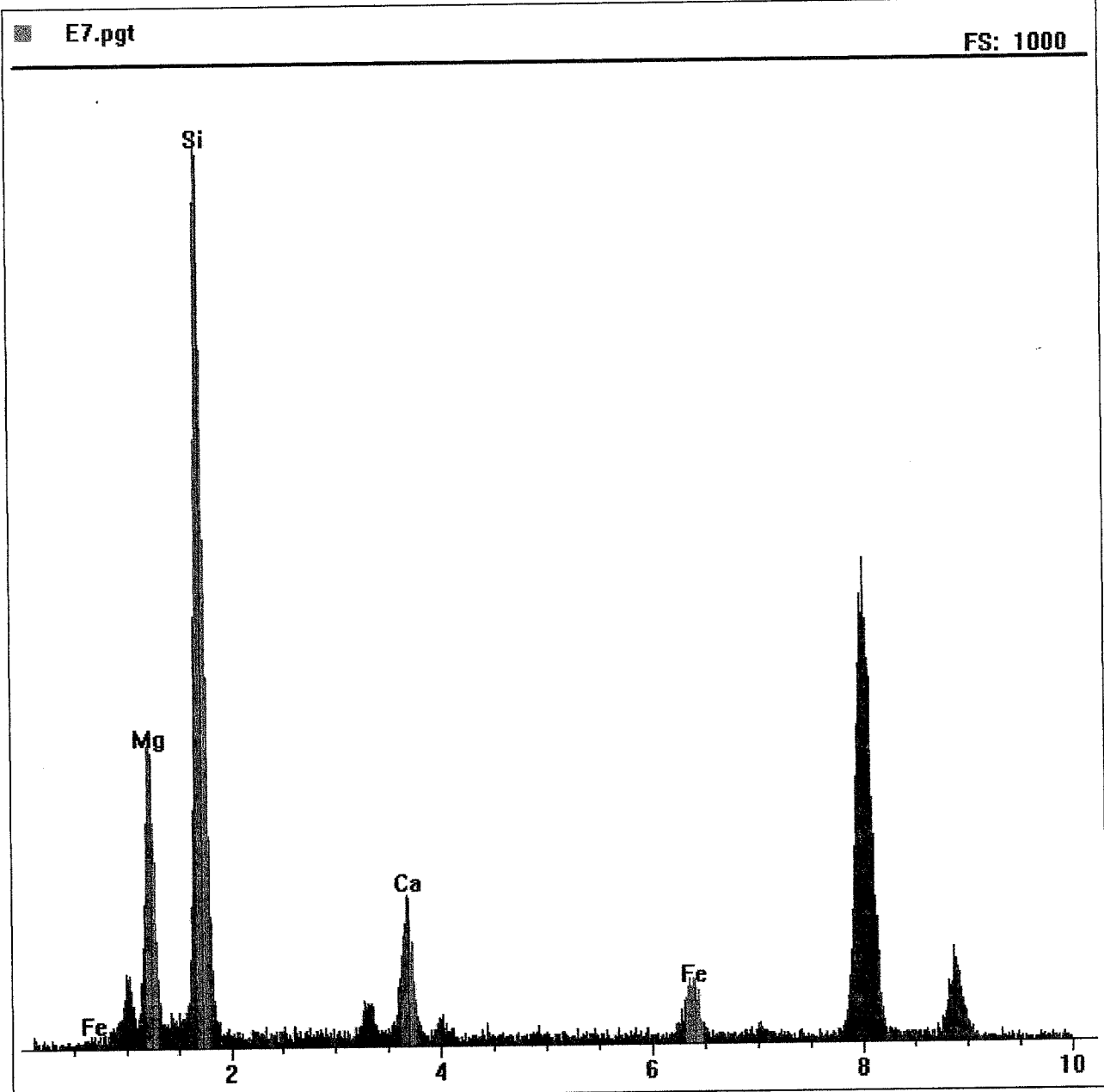

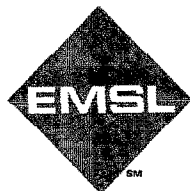

Energy Dispersive X-Ray Analysis

Qualitative Spectrum


EMSL ANALYTICAL, INC.


File: F:\Documen...records\PGT Files\EMSL27-2\EMSL27-2 2013\SW21281\E8.pgt
Collected: September 24, 2013 10:53:11

Live Time: 170.70
Beam Voltage: 20.00

Count Rate: 1239
Beam Current: 2.00

Dead Time: 15.00 %
Takeoff Angle: 60.98

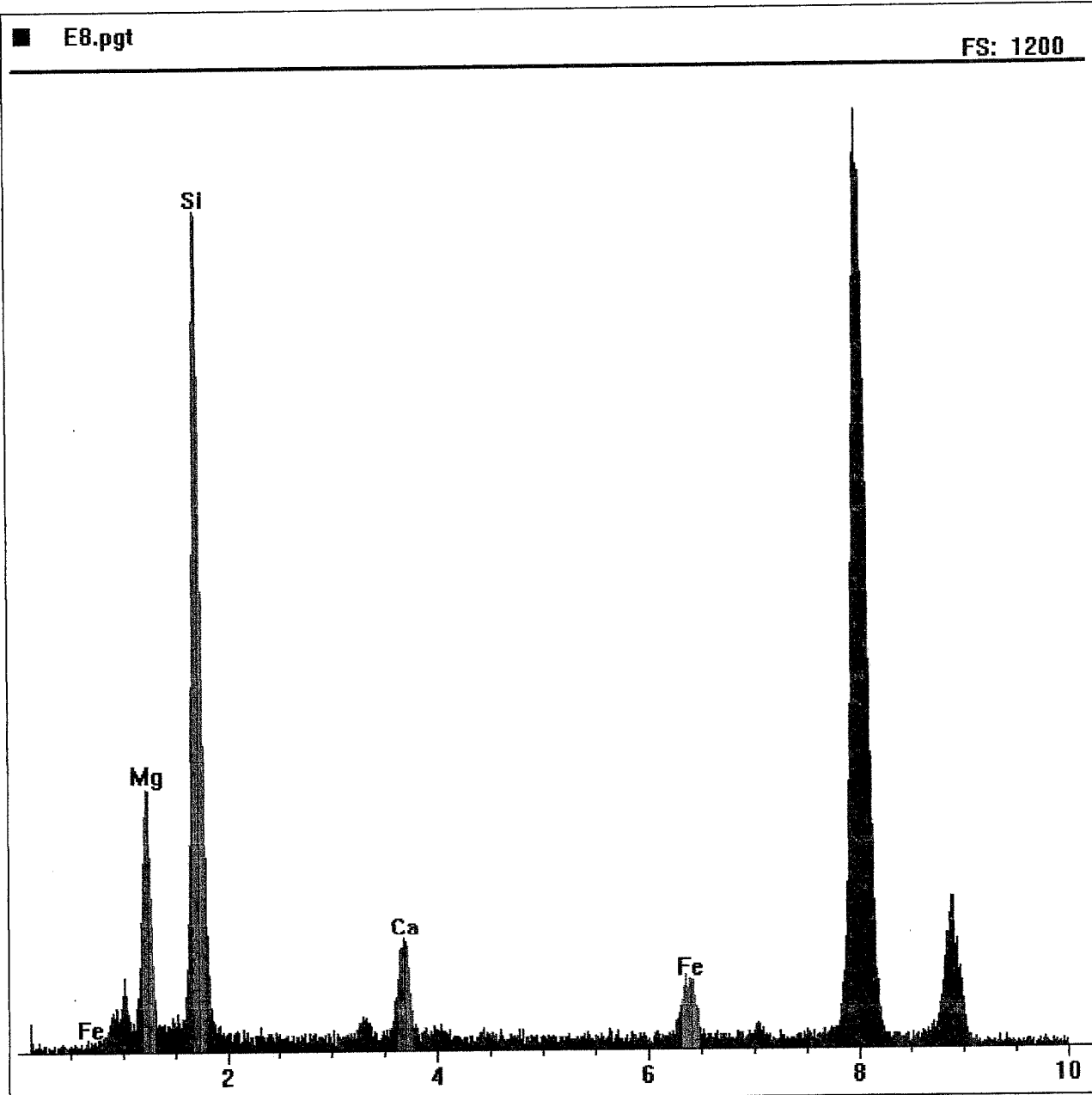

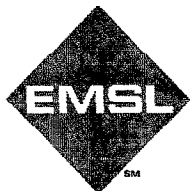

Energy Dispersive X-Ray Analysis Qualitative Spectrum


EMSL ANALYTICAL, INC.


File: F:\Documen...ecords\PGT Files\EMSL27-2\EMSL27-2 2013\SW21281\E9.pgt
Collected: September 24, 2013 10:53:11

Live Time: 228.41
Beam Voltage: 20.00

Count Rate: 742
Beam Current: 2.00

Dead Time: 11.62 %
Takeoff Angle: 60.98

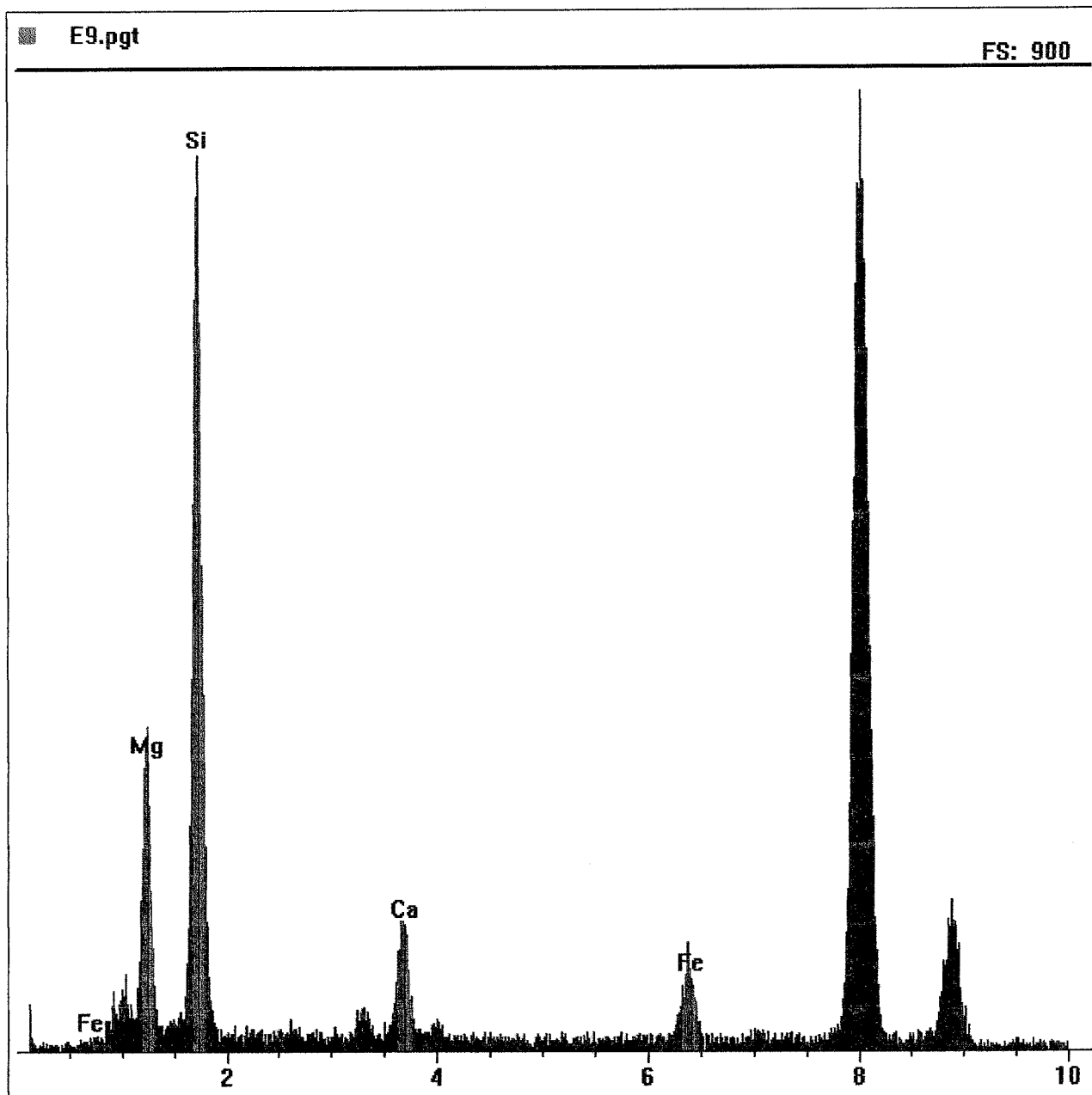

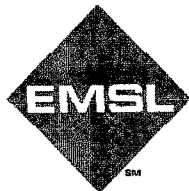

Energy Dispersive X-Ray Analysis

Qualitative Spectrum


EMSL ANALYTICAL, INC.


File: F:\Documen...ords\PGT Files\EMSL27-2\EMSL27-2 2013\SW21281\E10.pgt
Collected: September 24, 2013 10:53:11

Live Time: 122.24
Beam Voltage: 20.00

Count Rate: 1532
Beam Current: 2.00

Dead Time: 17.59 %
Takeoff Angle: 60.98

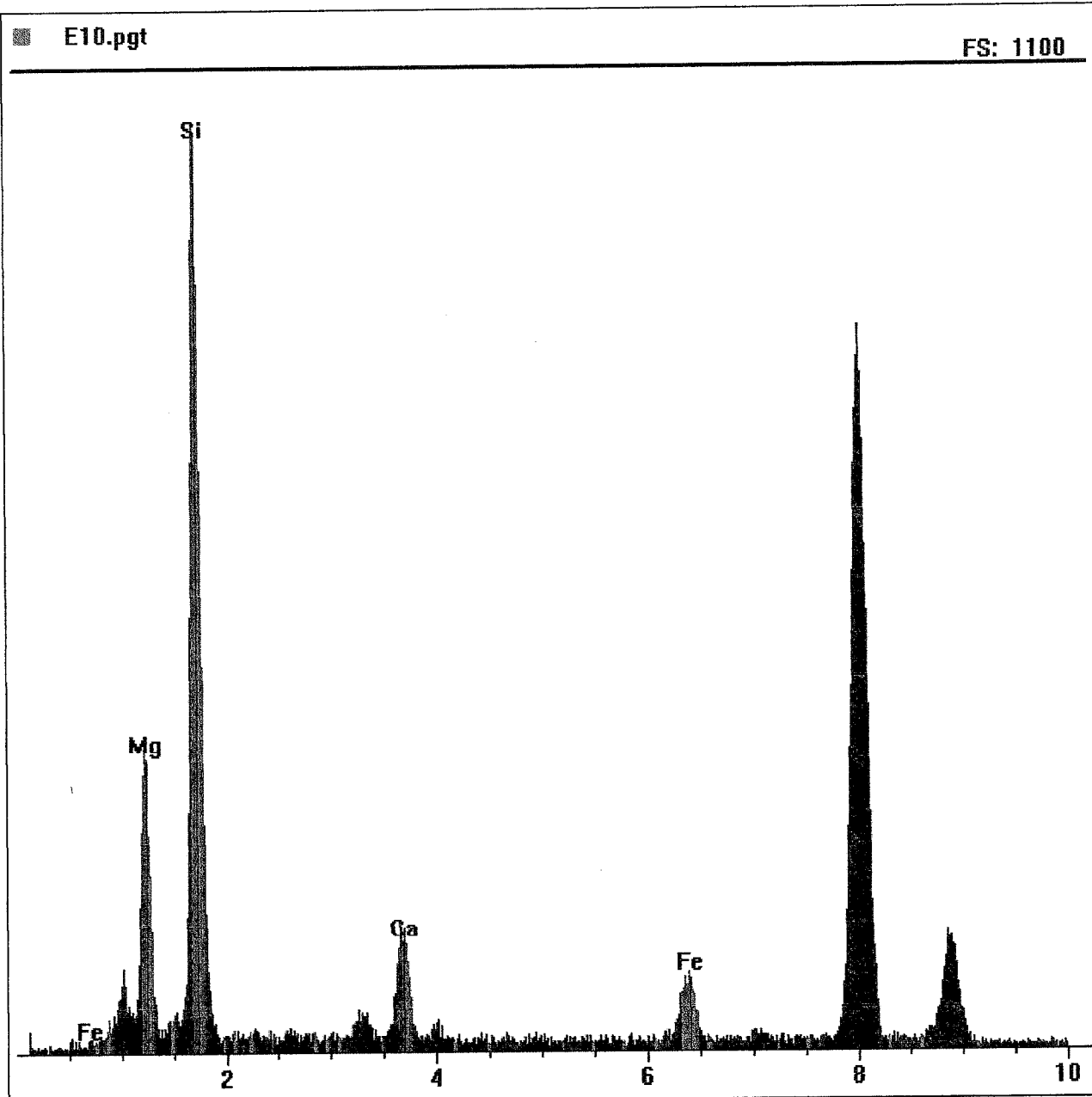

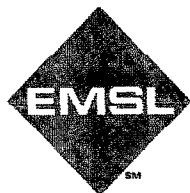

Energy Dispersive X-Ray Analysis

Qualitative Spectrum


EMSL ANALYTICAL, INC.


File: F:\Documen...ecords\PGT Files\EMSL27-2\EMSL27-2 2013\SW01TR1\G1.pgt
Collected: September 24, 2013 11:11:27

Live Time: 21.99
Beam Voltage: 20.00

Count Rate: 5806
Beam Current: 2.00

Dead Time: 48.53 %
Takeoff Angle: 60.98

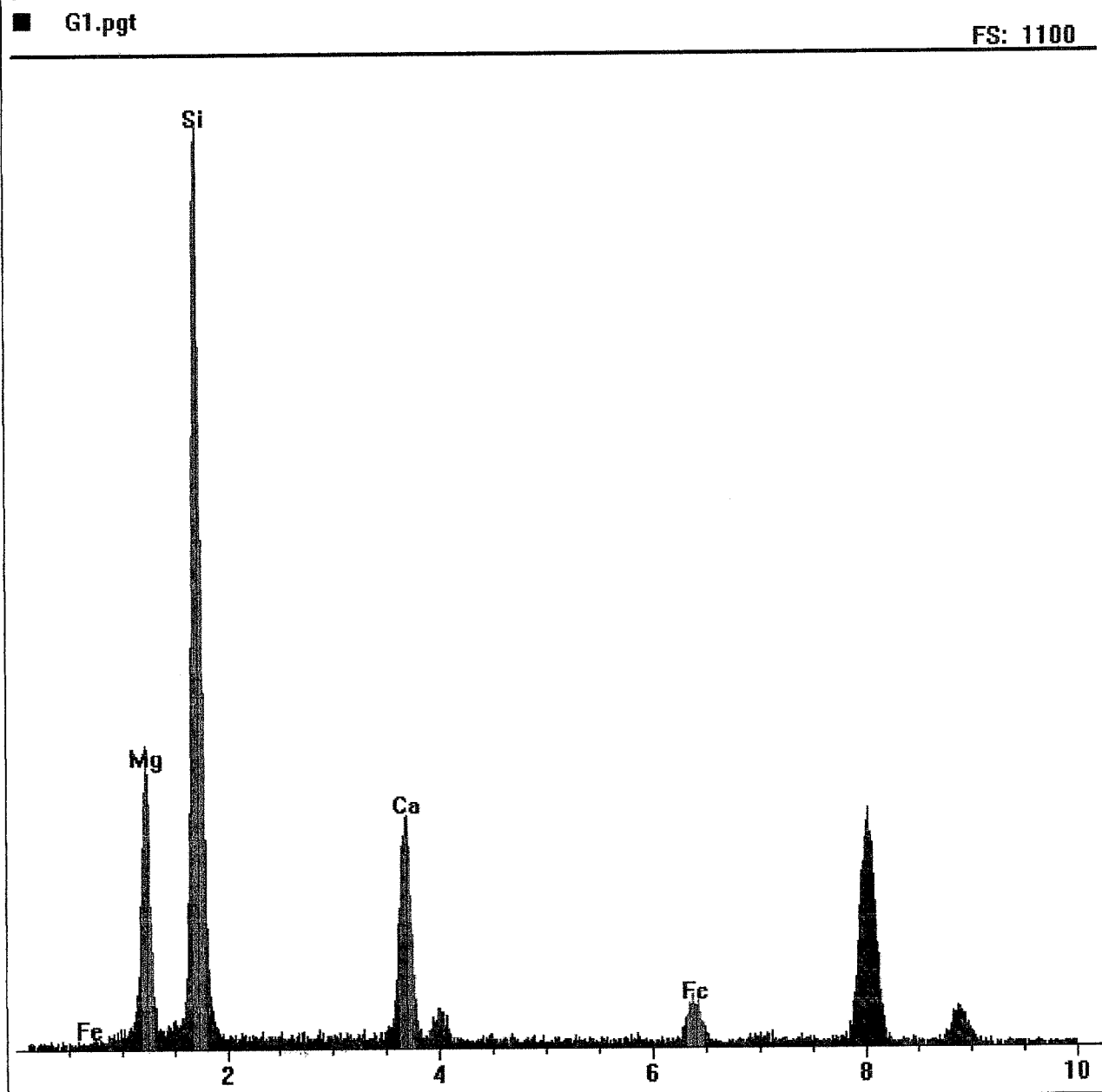

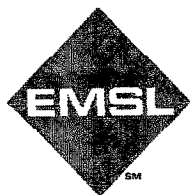

Energy Dispersive X-Ray Analysis

Qualitative Spectrum


EMSL ANALYTICAL, INC.


File: F:\Documen...records\PGT Files\EMSL27-2\EMSL27-2 2013\SW01TR1\G2.pgt
Collected: September 24, 2013 11:11:27

Live Time: 104.25
Beam Voltage: 20.00

Count Rate: 1461
Beam Current: 2.00

Dead Time: 16.74 %
Takeoff Angle: 60.98

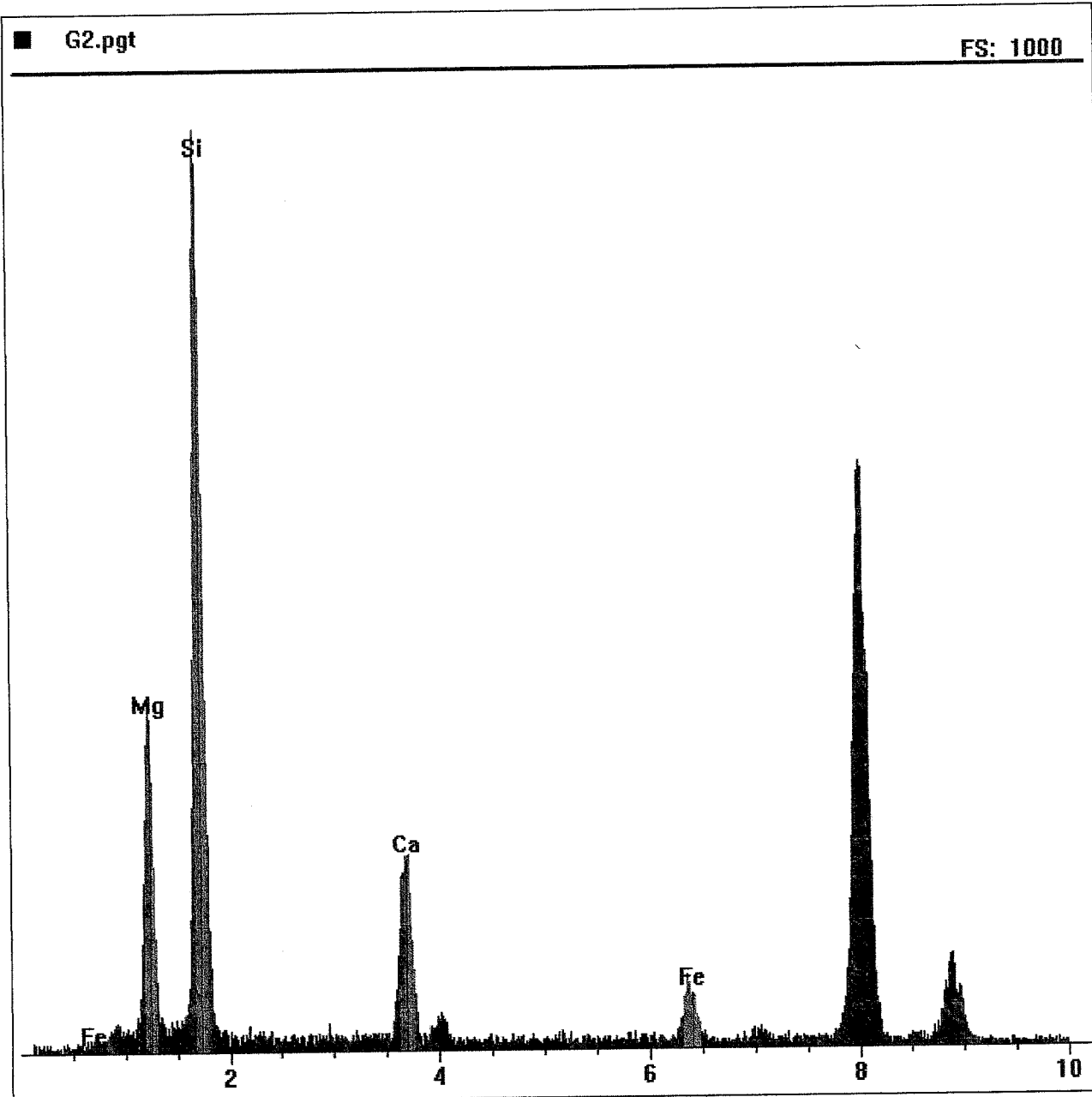

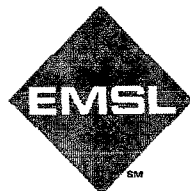

Energy Dispersive X-Ray Analysis

Qualitative Spectrum

EMSL ANALYTICAL, INC.

File: F:\Documen...ecords\PGT Files\EMSL27-2\EMSL27-2 2013\SW01TR1\G3.pgt

Collected: September 24, 2013 11:11:27

Live Time: 45.92


Count Rate: 3193


Dead Time: 30.84 %

Beam Voltage: 20.00

Beam Current: 2.00

Takeoff Angle: 60.98

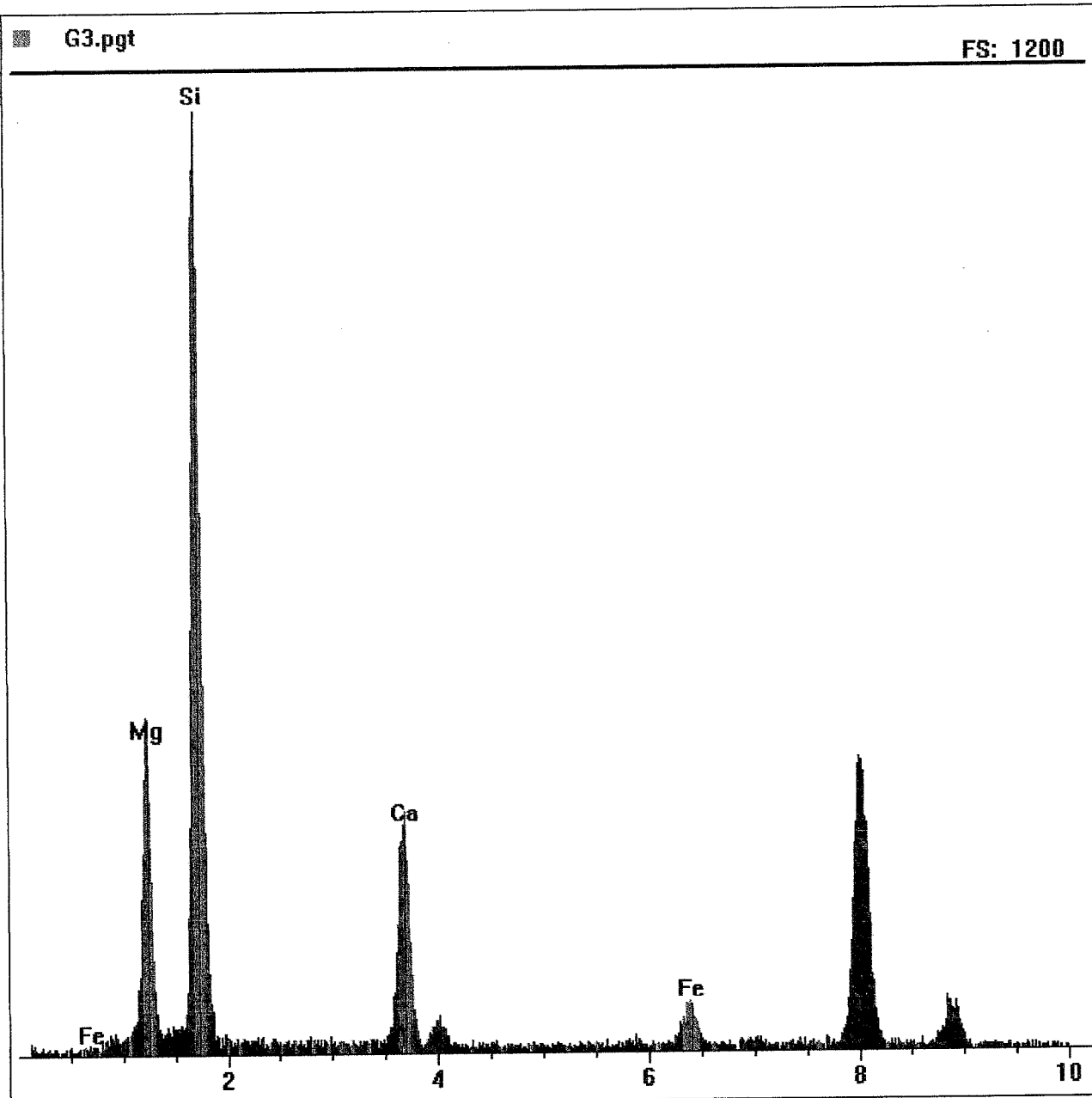

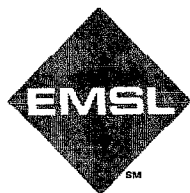

Energy Dispersive X-Ray Analysis

Qualitative Spectrum


EMSL ANALYTICAL, INC.


File: F:\Documen...ecords\PGT Files\EMSL27-2\EMSL27-2 2013\SW01TR1\G4.pgt
Collected: September 24, 2013 11:17:49

Live Time: 65.85
Beam Voltage: 20.00

Count Rate: 2492
Beam Current: 2.00

Dead Time: 25.48 %
Takeoff Angle: 60.98

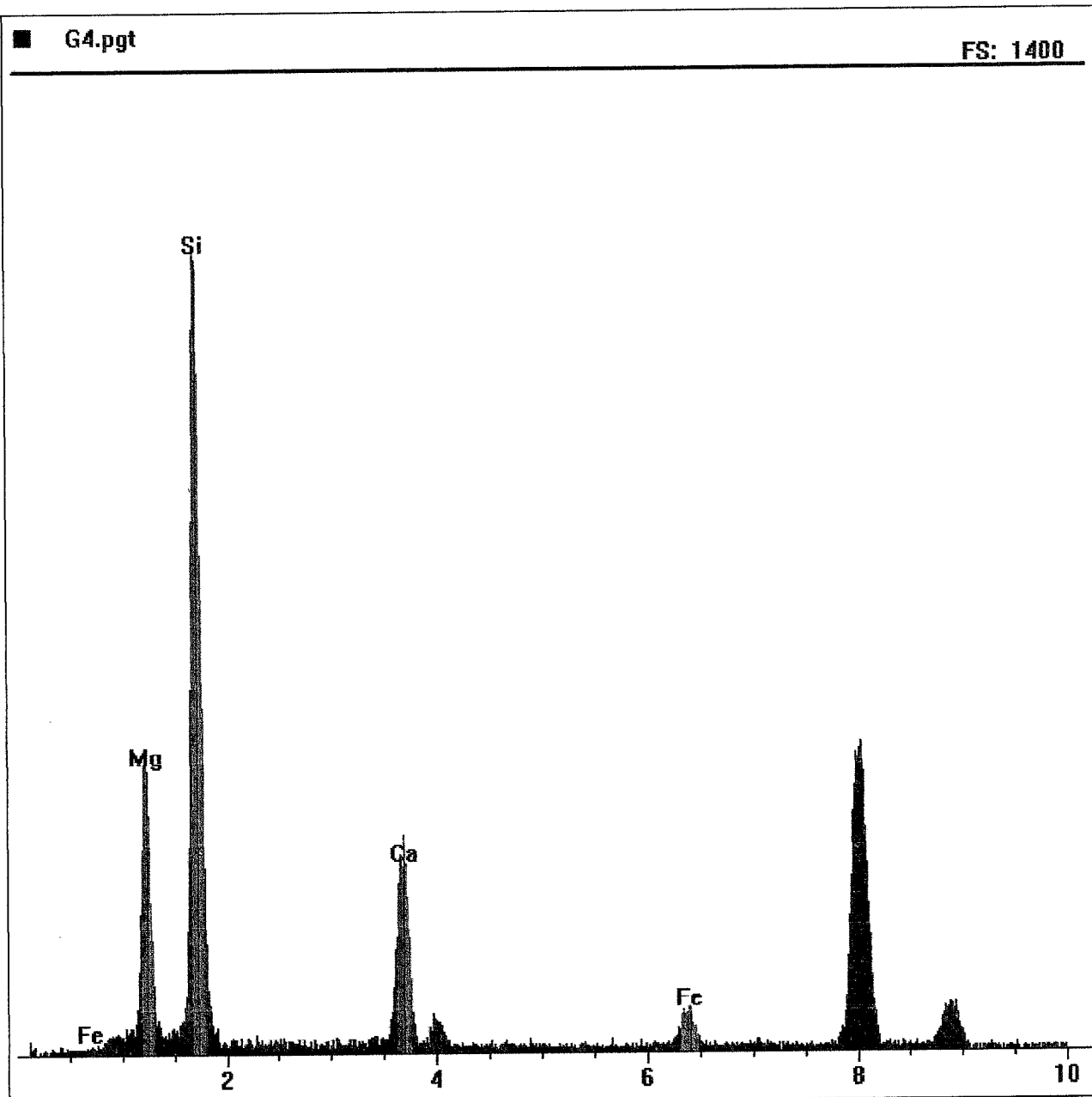

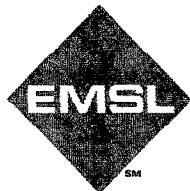

Energy Dispersive X-Ray Analysis

Qualitative Spectrum


EMSL ANALYTICAL, INC.


File: F:\Documen...ecords\PGT Files\EMSL27-2\EMSL27-2 2013\SW01TR1\G5.pgt
Collected: September 24, 2013 11:17:49

Live Time: 65.85
Beam Voltage: 20.00

Count Rate: 2492
Beam Current: 2.00

Dead Time: 25.48 %
Takeoff Angle: 60.98

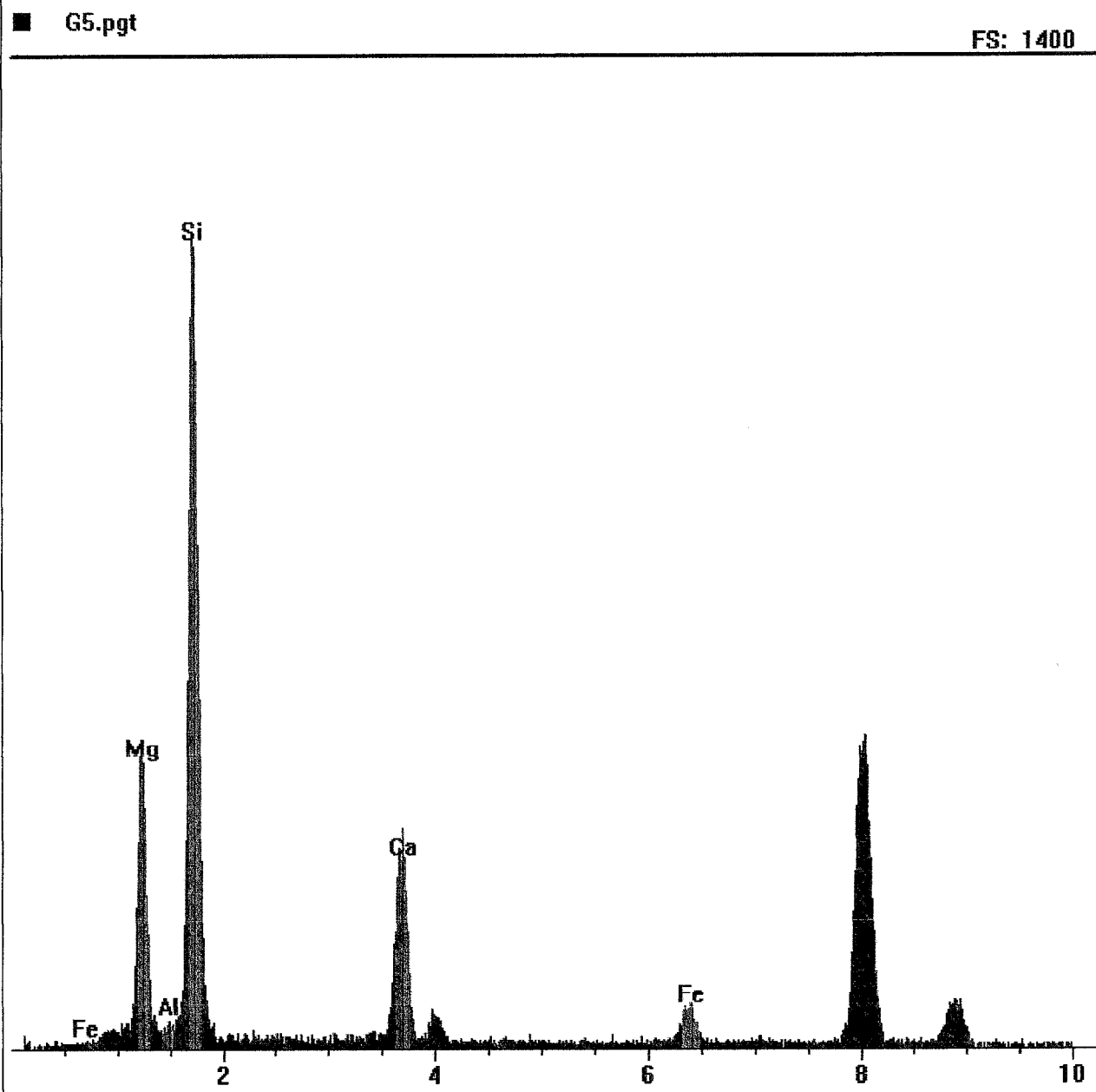

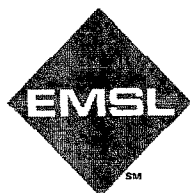

Energy Dispersive X-Ray Analysis

Qualitative Spectrum

EMSL ANALYTICAL, INC.

File: F:\Documen...ecords\PGT Files\EMSL27-2\EMSL27-2 2013\SW01TR1\G6.pgt

Collected: September 24, 2013 11:17:49

Live Time: 43.25


Count Rate: 2735

Dead Time: 27.38 %

Beam Voltage: 20.00

Beam Current: 2.00

Takeoff Angle: 60.98

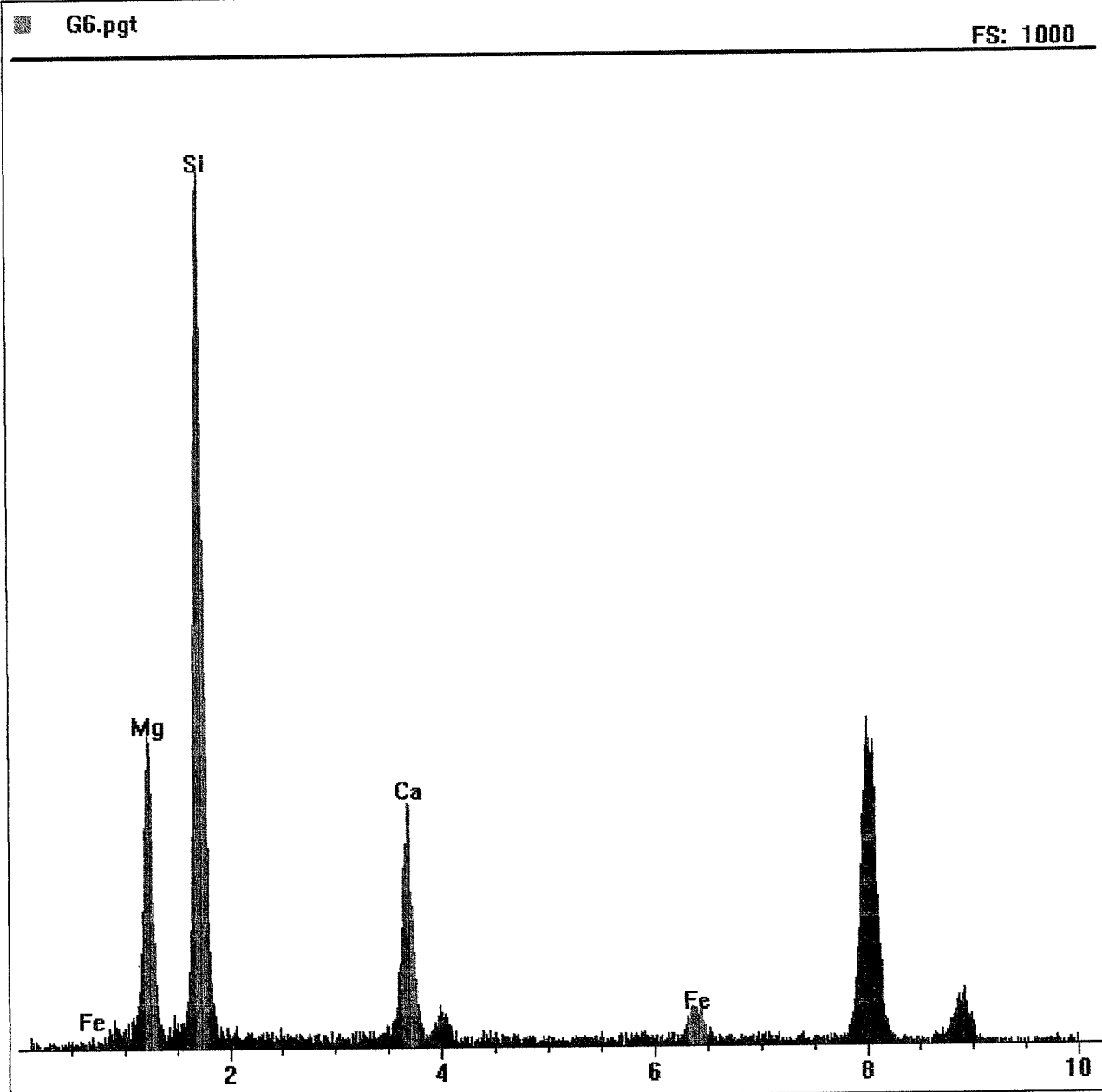



Energy Dispersive X-Ray Analysis Qualitative Spectrum


EMSL ANALYTICAL, INC.

File: F:\Documen...ecords\PGT Files\EMSL27-2\EMSL27-2 2013\SW01TR1\G7.pgt
Collected: September 24, 2013 11:21:56

Live Time: 104.15
Beam Voltage: 20.00

Count Rate: 1460
Beam Current: 2.00

Dead Time: 16.76 %
Takeoff Angle: 60.98

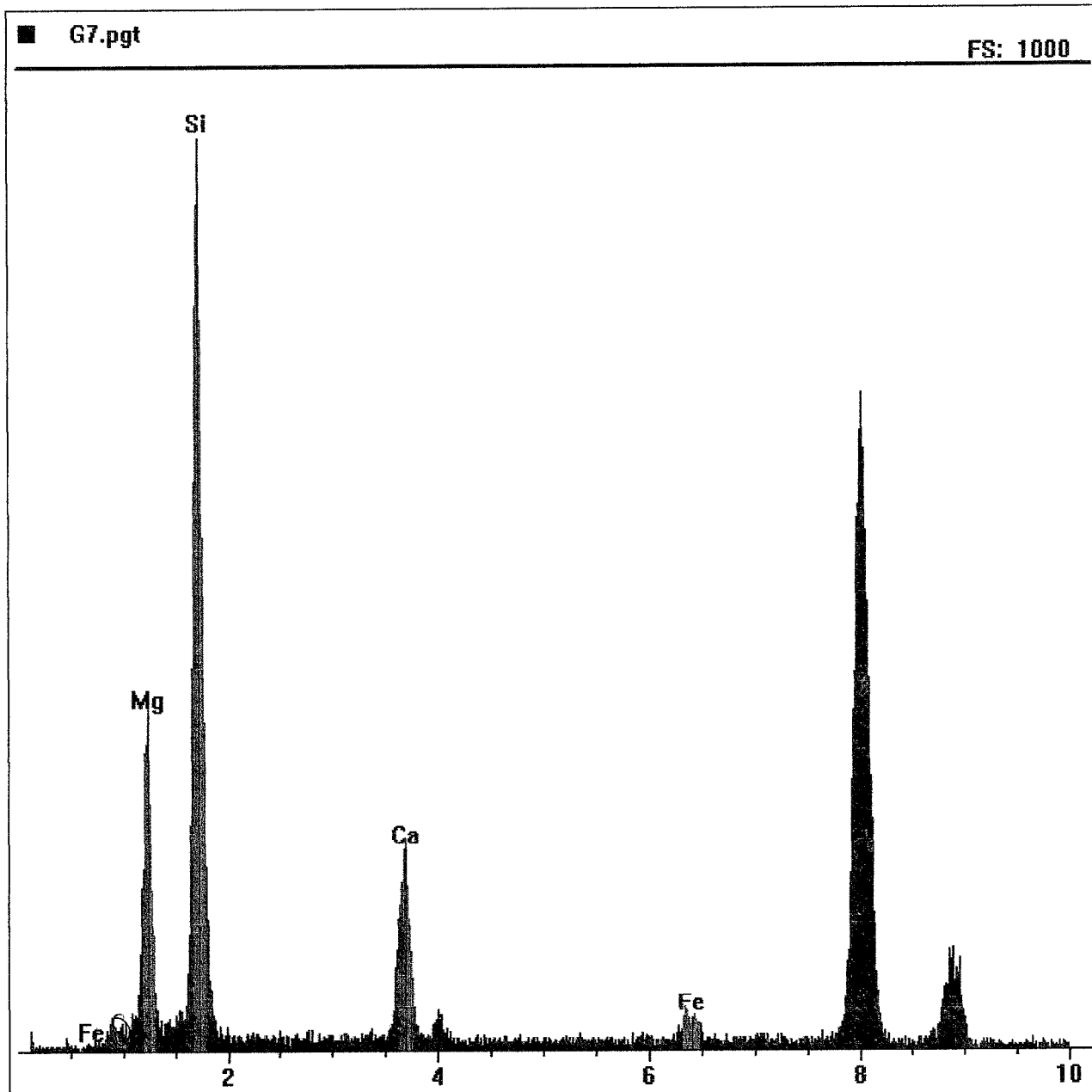



Energy Dispersive X-Ray Analysis

Qualitative Spectrum

EMSL ANALYTICAL, INC.

File: F:\Documen...ecords\PGT Files\EMSL27-2\EMSL27-2 2013\SW01TR1\G8.pgt

Collected: September 24, 2013 11:21:56

Live Time: 48.97


Count Rate: 4057


Dead Time: 37.03 %

Beam Voltage: 20.00

Beam Current: 2.00

Takeoff Angle: 60.98

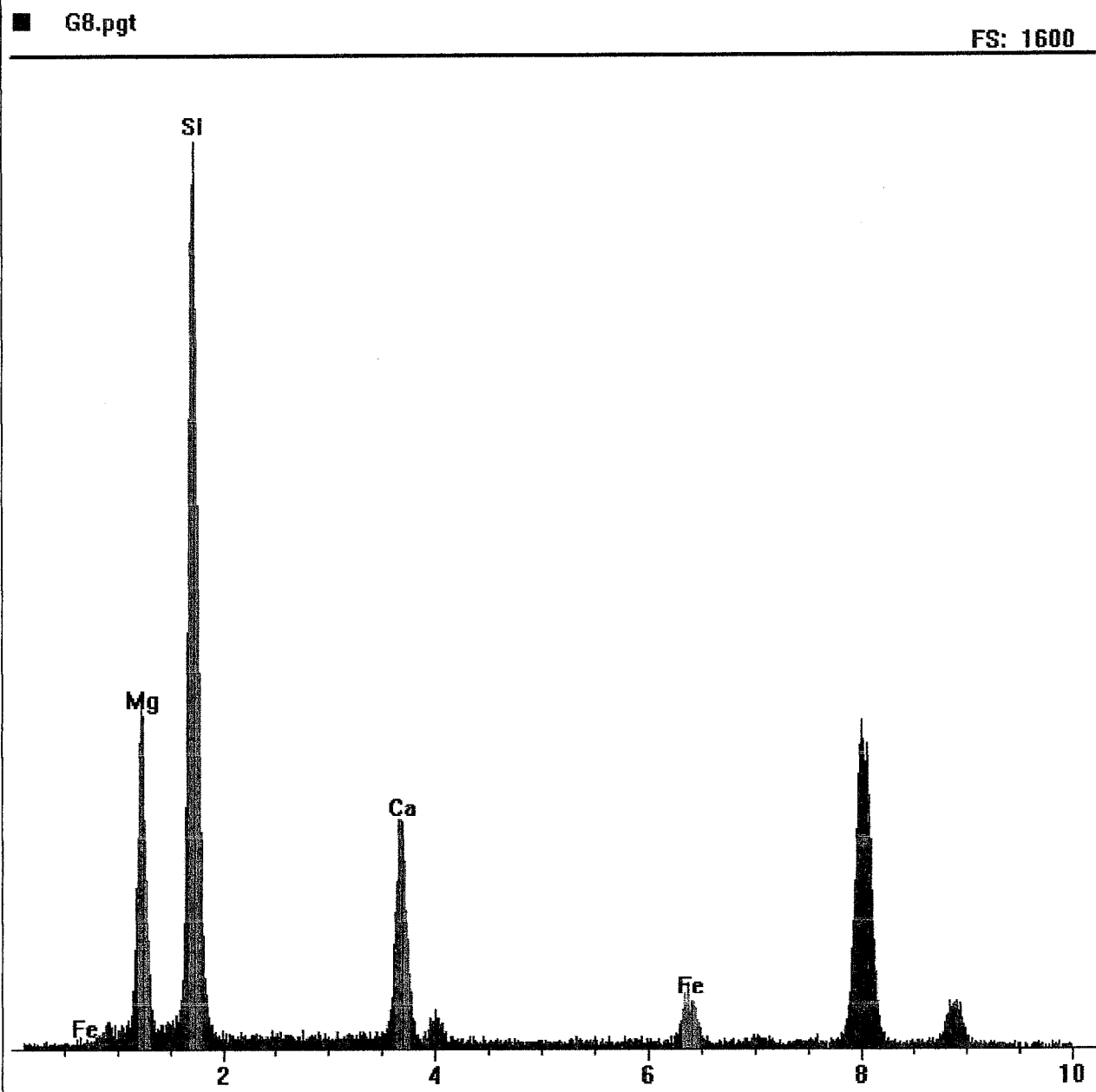

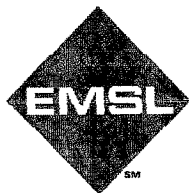

Energy Dispersive X-Ray Analysis

Qualitative Spectrum

EMSL ANALYTICAL, INC.

File: F:\Documen...records\PGT Files\EMSL27-2\EMSL27-2 2013\SW01TR1\G9.pgt

Collected: September 24, 2013 11:21:56

Live Time: 34.59


Count Rate: 3998


Dead Time: 36.75 %

Beam Voltage: 20.00

Beam Current: 2.00

Takeoff Angle: 60.98

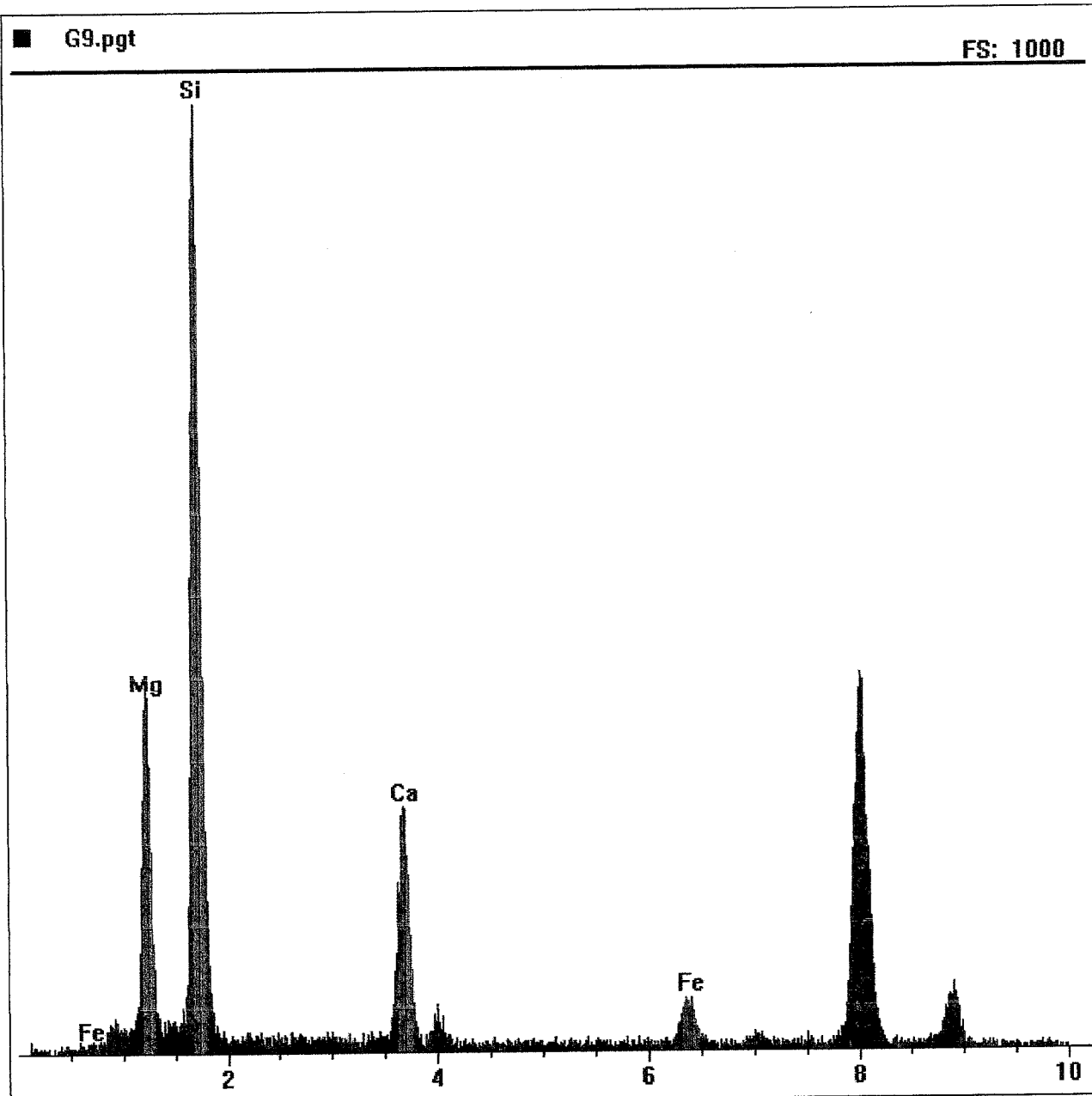

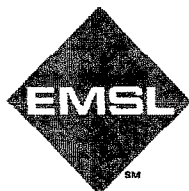

Energy Dispersive X-Ray Analysis

Qualitative Spectrum

EMSL ANALYTICAL, INC.

File: F:\Documen...ords\PGT Files\EMSL27-2\EMSL27-2 2013\SW01TR1\G10.pgt
Collected: September 24, 2013 11:21:56

Live Time: 12.47

Count Rate: 11299

Dead Time: 72.57 %


Beam Voltage: 20.00

Beam Current: 2.00

Takeoff Angle: 60.98

■ G10.pgt

FS: 1100

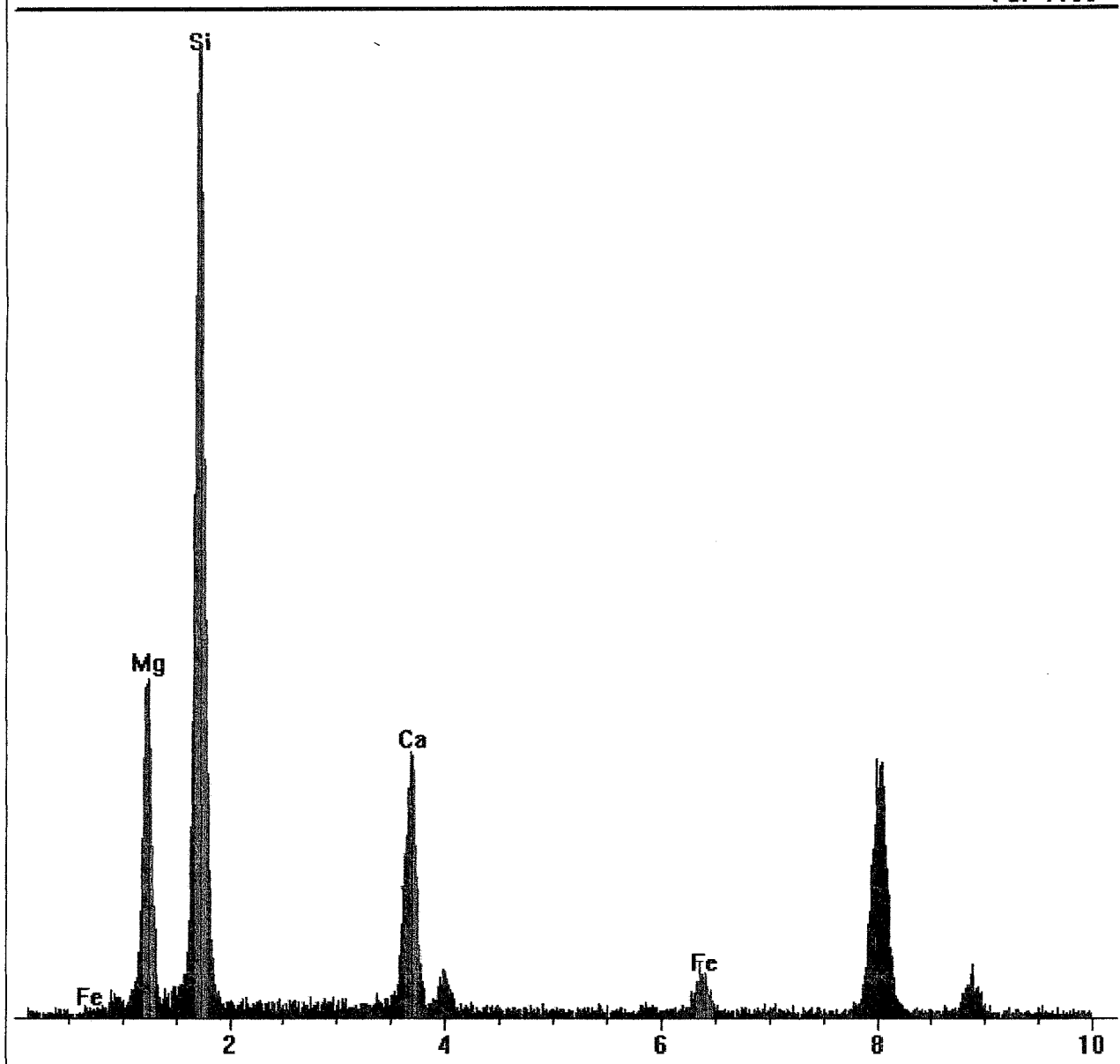



Energy Dispersive X-Ray Analysis Qualitative Spectrum


EMSL ANALYTICAL, INC.

File: F:\Documen...2\EMSL27-2 2013\SRM1867 Actinolite\SRM1867_Act_05t.pgt
Collected: September 27, 2013 08:33:01

Live Time: 208.79
Beam Voltage: 20.00

Count Rate: 1240
Beam Current: 2.00

Dead Time: 14.68 %
Takeoff Angle: 60.98

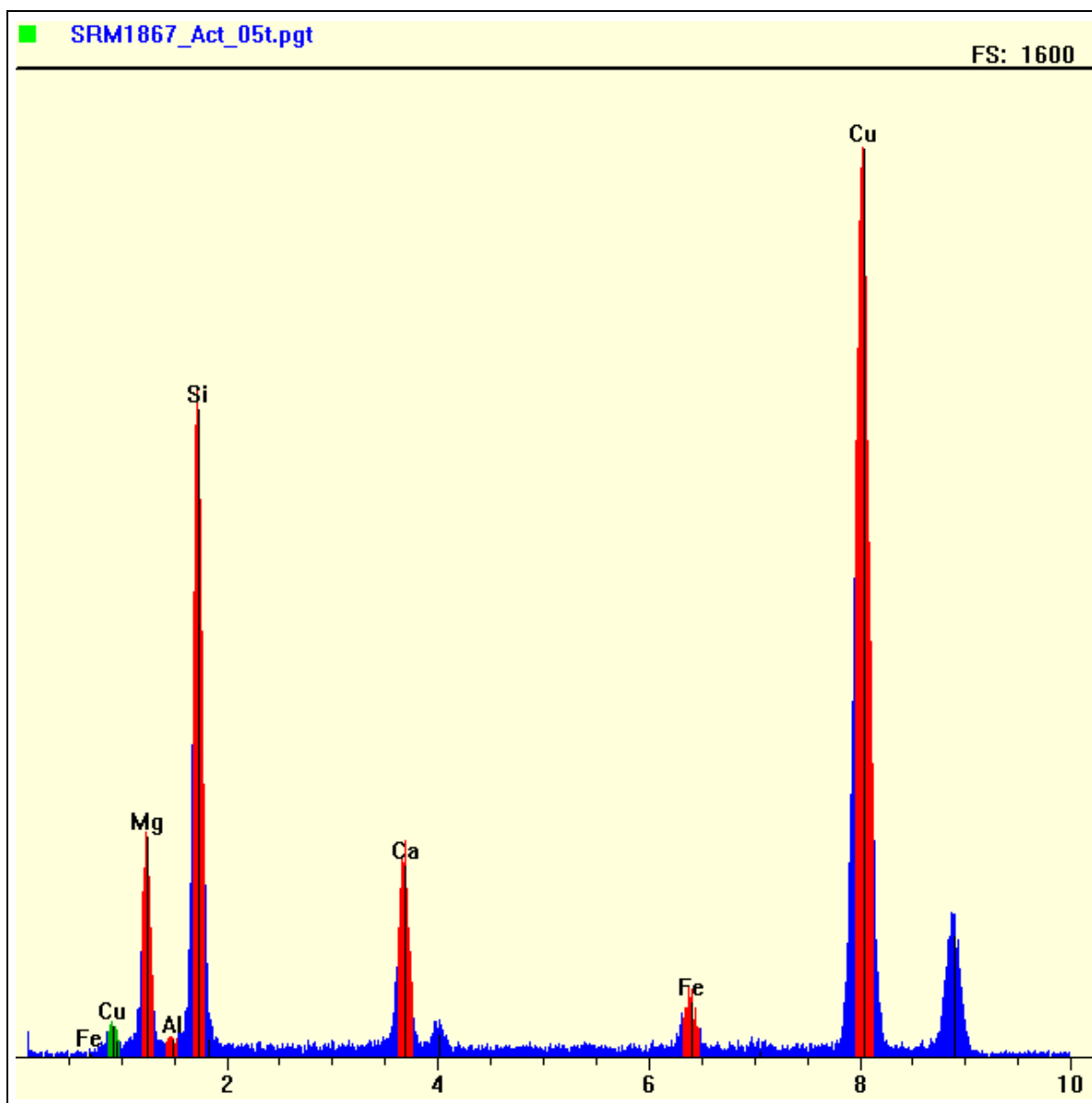

EDXA of NIST Actinolite Image 29019



Energy Dispersive X-Ray Analysis Qualitative Spectrum


EMSL ANALYTICAL, INC.

File: F:\Documen...-2\EMSL27-2 2013\SRM1867 Actinolite\SRM1867_Act_15.pgt
Collected: September 27, 2013 08:33:01

Live Time: 408.63
Beam Voltage: 20.00

Count Rate: 777
Beam Current: 2.00

Dead Time: 10.72 %
Takeoff Angle: 60.98

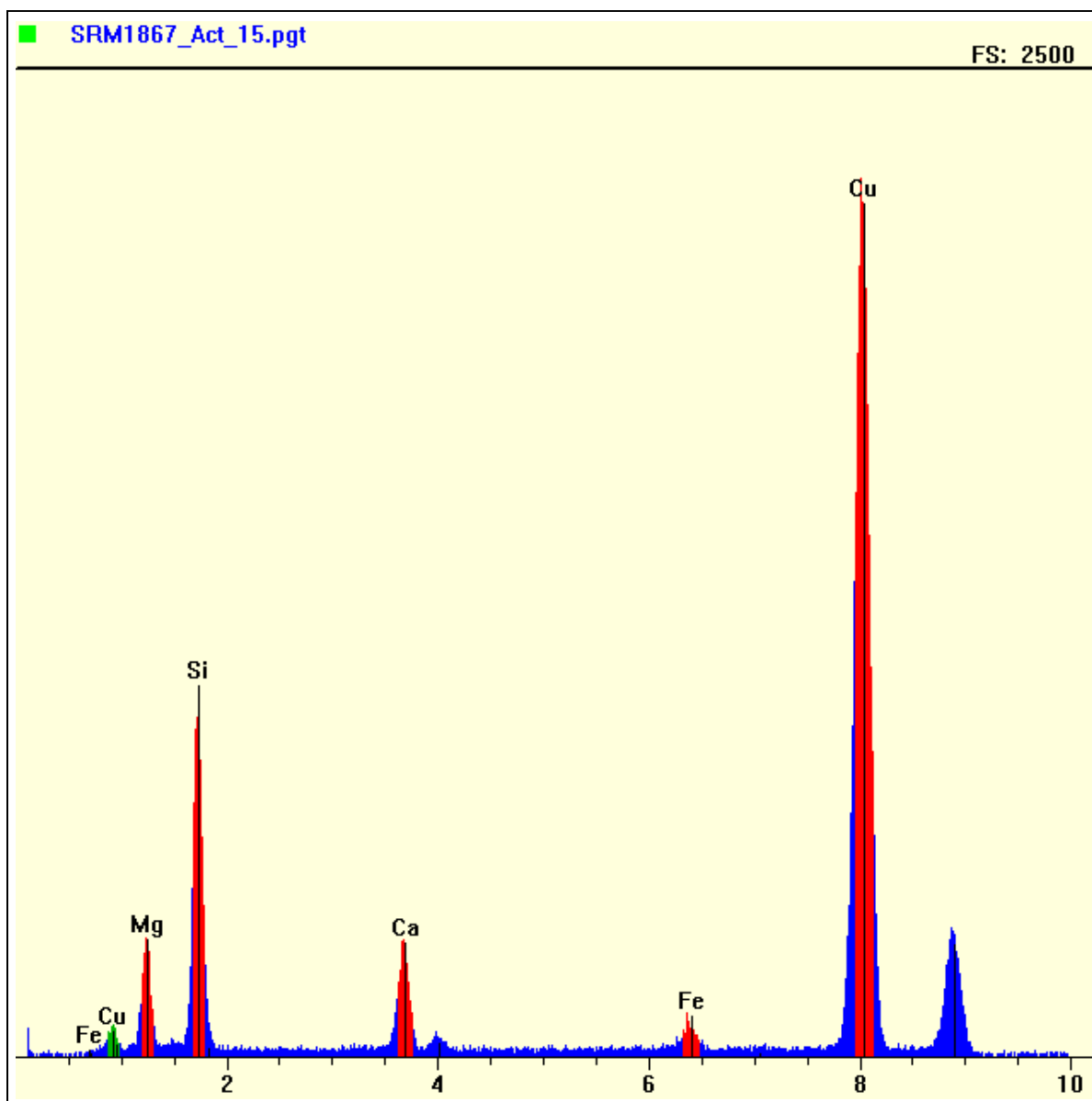

EDXA of NIST Actinolite Image 29031

EMSL Analytical Inc. 107 West 4th Street, Libby, MT, 59923
TEL: (406) 293-9066 • FAX: (406) 293-7016 • www.EMSL.com
Locally Focused ... Nationally Recognized



Energy Dispersive X-Ray Analysis Qualitative Spectrum


EMSL ANALYTICAL, INC.

File: F:\Documen...-2\EMSL27-2 2013\SRM1867 Actinolite\SRM1867_Act_19.pgt
Collected: September 27, 2013 08:33:01

Live Time: 104.32
Beam Voltage: 20.00

Count Rate: 1895
Beam Current: 2.00

Dead Time: 20.46 %
Takeoff Angle: 60.98

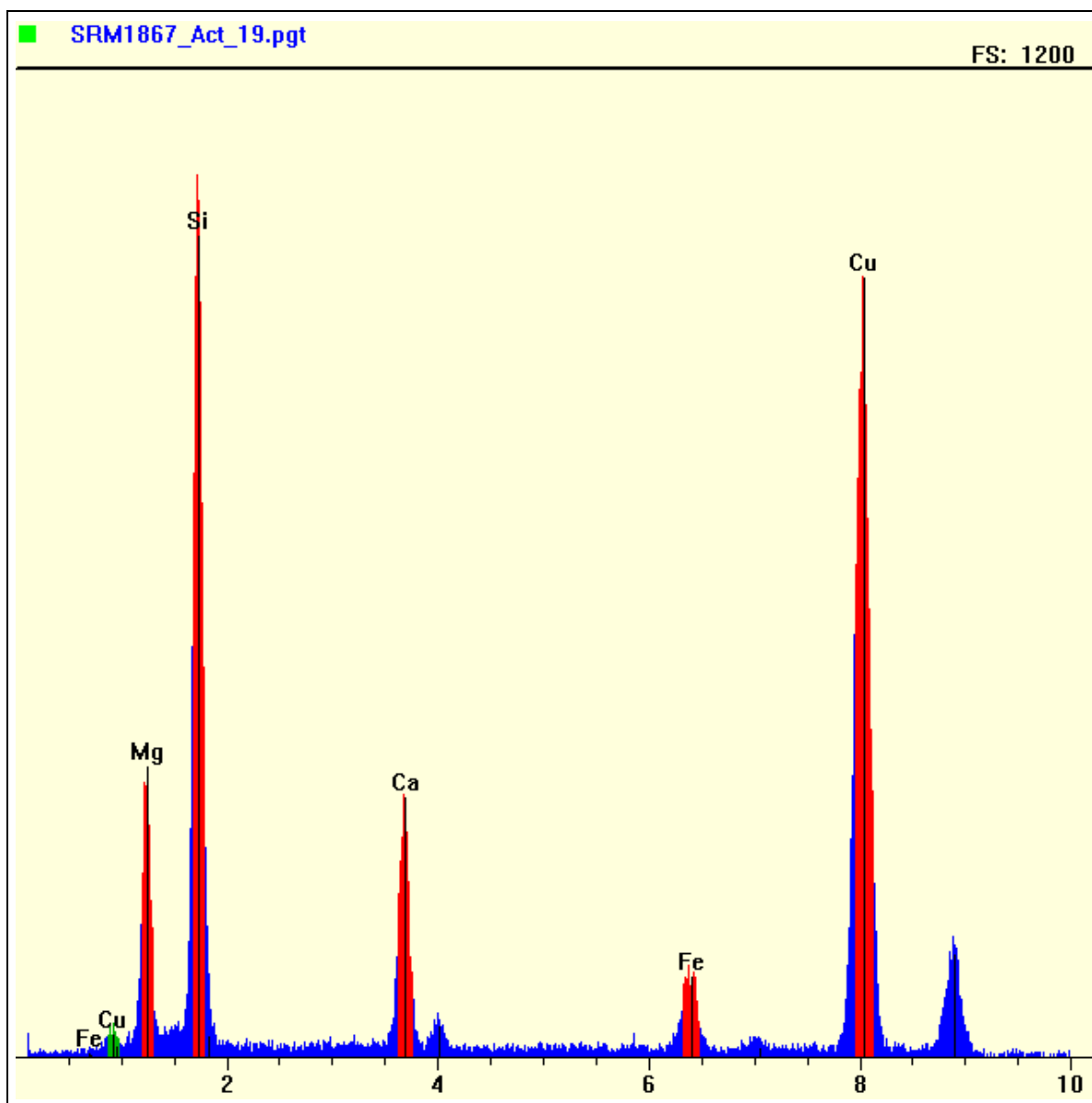

EDXA of NIST Actinolite Image 29035



Energy Dispersive X-Ray Analysis Qualitative Spectrum


EMSL ANALYTICAL, INC.

File: F:\Documen...L27-2\EMSL27-2 2013\Diopside_01\Diopside_093013_05.pgt
Collected: September 30, 2013 10:16:44

Live Time: 108.71
Beam Voltage: 20.00

Count Rate: 1426
Beam Current: 2.00

Dead Time: 16.43 %
Takeoff Angle: 60.98

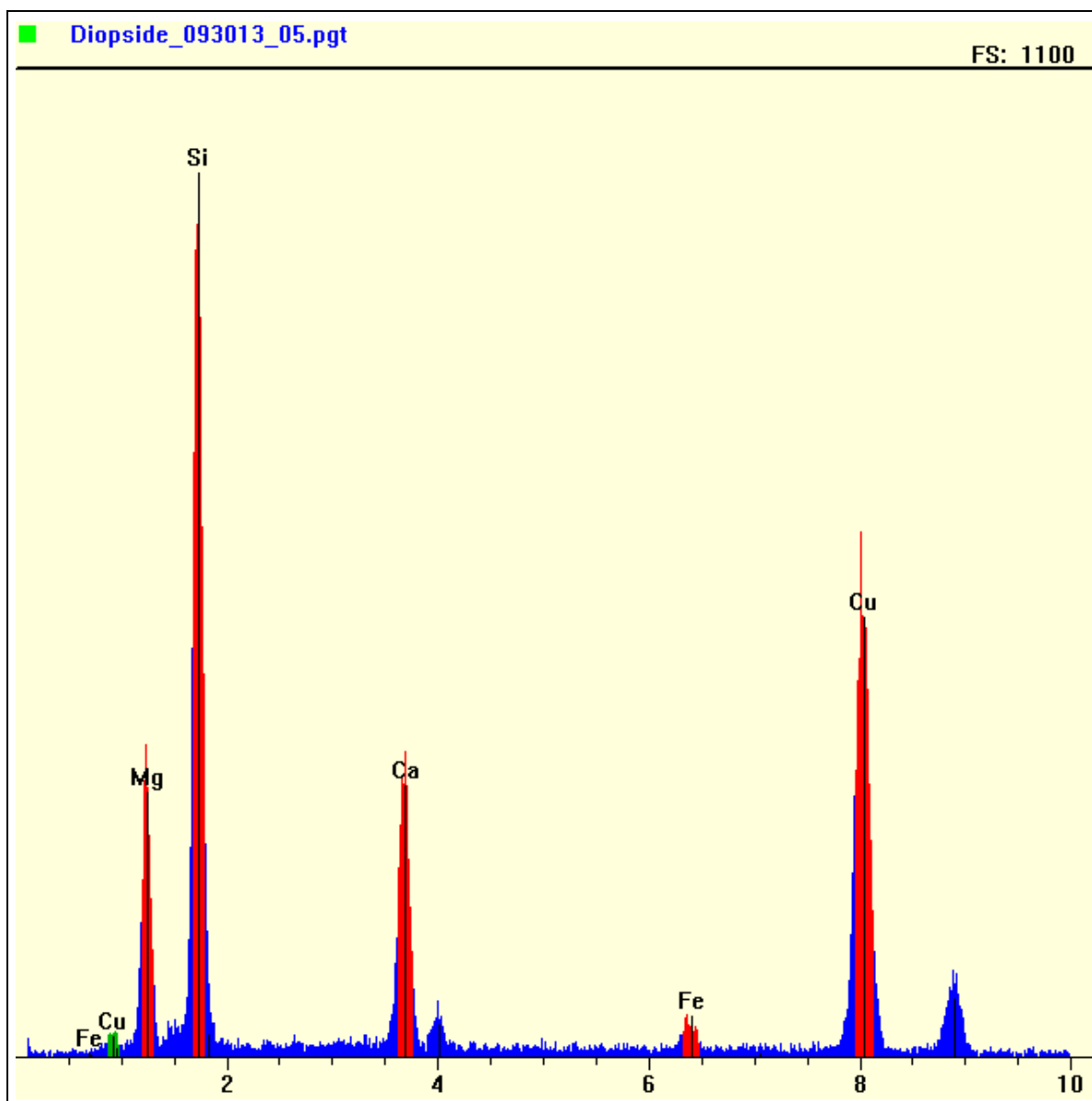

EDXA of Diopside Image 29041



Energy Dispersive X-Ray Analysis Qualitative Spectrum


EMSL ANALYTICAL, INC.

File: F:\Documen...L27-2\EMSL27-2 2013\Diopside_01\Diopside_093013_10.pgt
Collected: September 30, 2013 10:16:44

Live Time: 150.73
Beam Voltage: 20.00

Count Rate: 1174
Beam Current: 2.00

Dead Time: 13.94 %
Takeoff Angle: 60.98

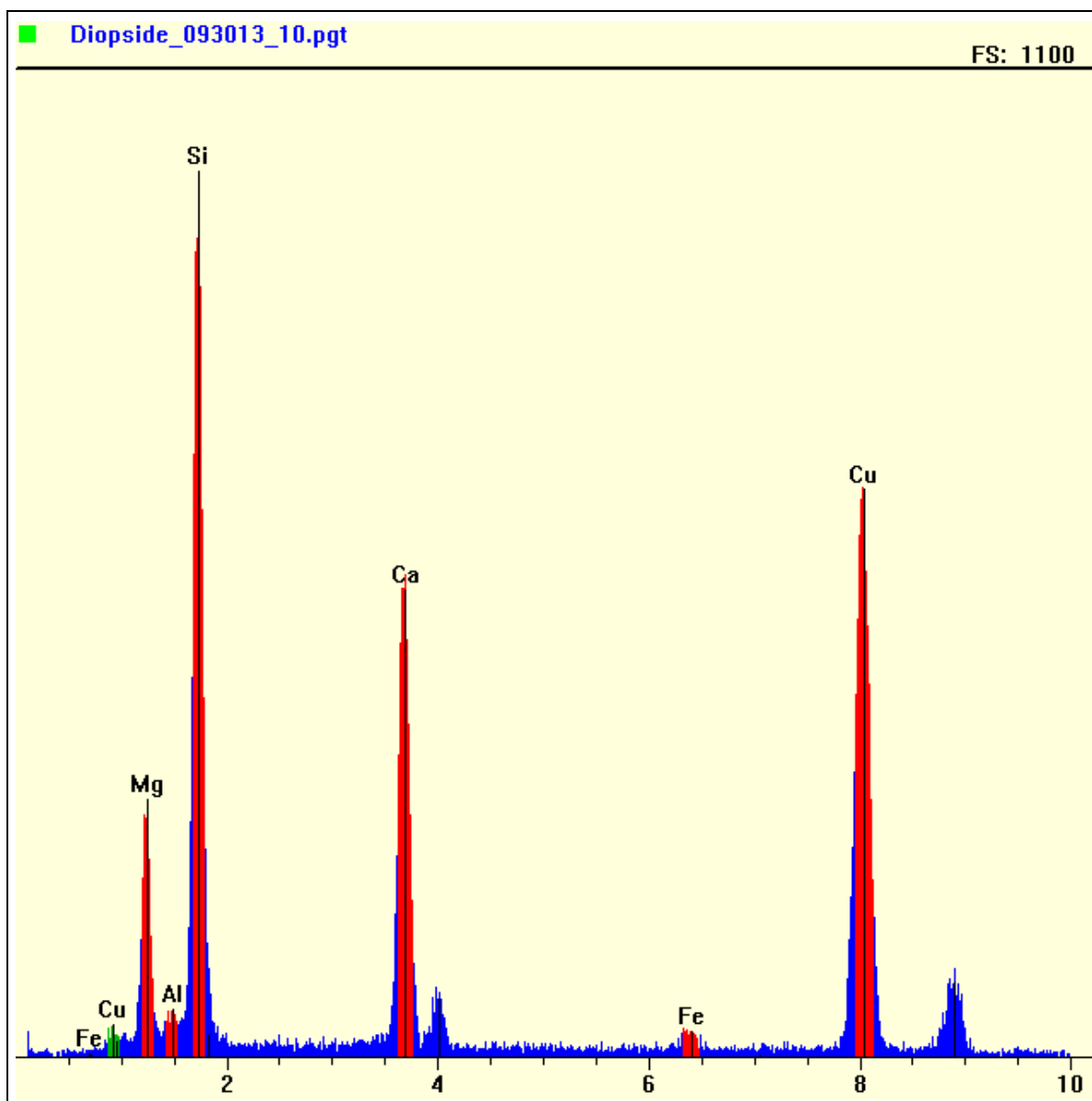

EDXA of Diopside Image 29047



Energy Dispersive X-Ray Analysis Qualitative Spectrum


EMSL ANALYTICAL, INC.

File: F:\Documen...L27-2\EMSL27-2 2013\Diopside_01\Diopside_093013_19.pgt
Collected: September 30, 2013 10:16:44

Live Time: 104.09
Beam Voltage: 20.00

Count Rate: 1515
Beam Current: 2.00

Dead Time: 17.75 %
Takeoff Angle: 60.98

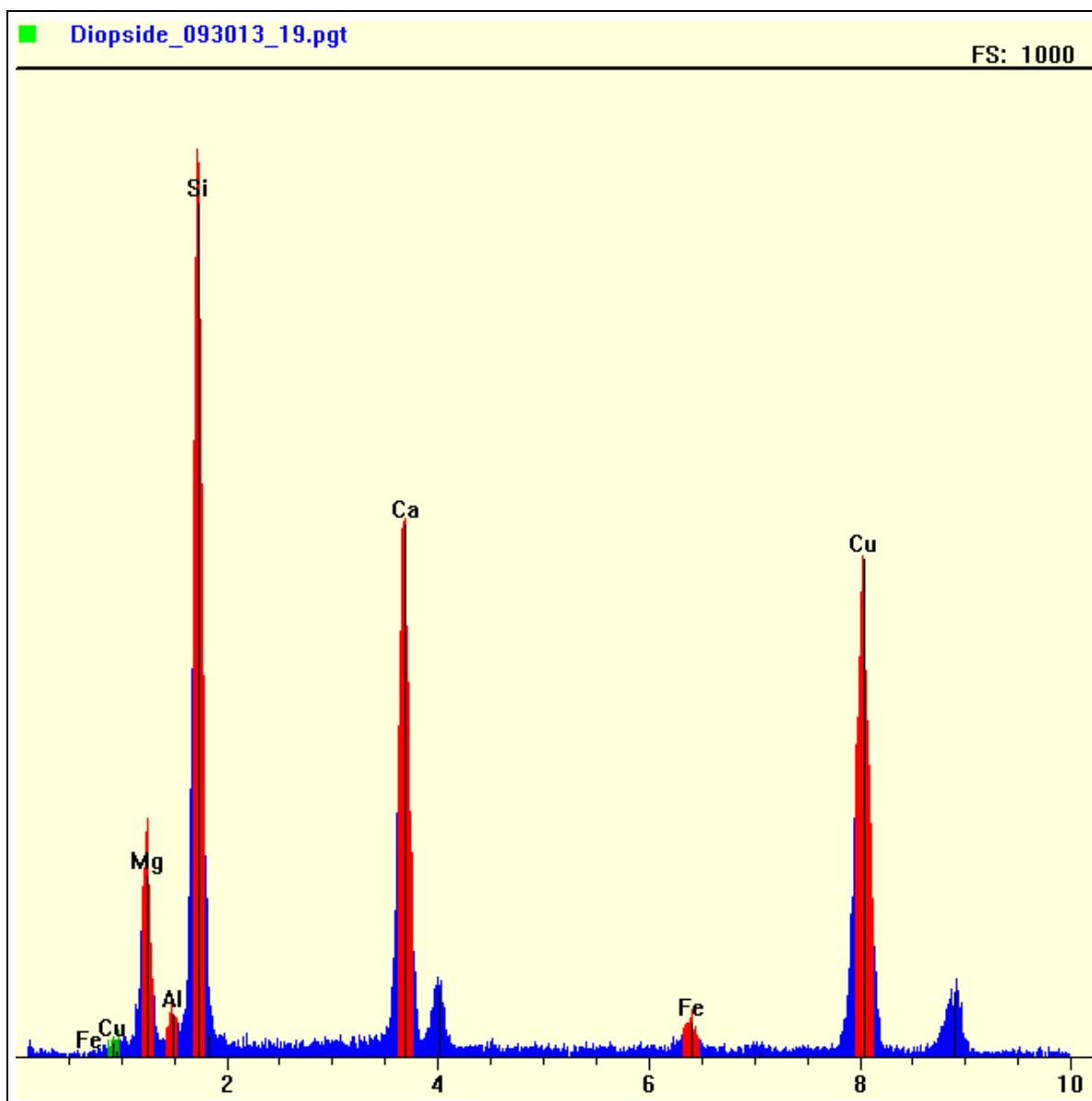

EDXA of Diopside Image 29057



AMPHIBOLE SAED INDEXING FORM

EMSL Order ID : <u>Actinolite</u>		DATE: <u>09/27/13</u>	
Indexing of negative number: <u>29019</u>		SCOPE # <u>27-2</u>	
Reference / Sample Number: <u>Actinolite</u>			
Preliminary ID: <u>Actinolite</u>		By: <u>EJWP</u>	
Using Camera Constant of: <u>22.82</u> mm Angstroms			
Determined from negative number: <u>29010</u>			
Measured Inter-row spacing:	<u>4.5</u>	mm	
Mean Distance between spots on Center row (d2):	<u>4.55</u>	mm	
Mean distance between spots on slant vector (d1):	<u>4.79</u>	mm	
	Calculated	Ref	- 5%
Inter-row spacing (Ångstroms)	5.07	5.300	5.035
d2 or hk0 (Camera K/zero row dist.)	5.02	5.099	4.844
d1 or hk1 (Camera K/slant vector dist.)	4.76	4.931	4.684
Ratio of hk0/hk1	1.05	1.03	0.979
Angle of Slant Vector (Measured)	68.0	68.0	64.600
	+ 5%		
			5.565
			5.354
			5.178
			1.082
			71.400
From SAED Reference Book, "unknown" diffraction pattern was found to be that of: <u>Actinolite (Reference JCPDS)</u> By: <u>EJWP</u> <i>and</i>			
With a Zone Axis of: [<u>3 1 2</u>]			
Preliminary Identification was: <input checked="" type="checkbox"/> CORRECT			
<input type="checkbox"/> INCORRECT			

percent accuracy to date: _____



AMPHIBOLE SAED INDEXING FORM

EMSL Order ID : Actinolite

DATE: 09/27/13

Indexing of negative number: 29031

SCOPE # 27-2

Reference / Sample Number: Actinolite

Preliminary ID: Actinolite

By: EJWP

Using Camera Constant of: 22.82 mm Angstroms

Determined from negative number: 29010

Measured Inter-row spacing:

4.25

mm

Mean Distance between spots on Center row (d2):

2.5

mm

Mean distance between spots on slant vector (d1):

5

mm

	Calculated	Ref	- 5%	+ 5%
Inter-row spacing (Ångstroms)	<u>5.37</u>	<u>5.300</u>	<u>5.035</u>	<u>5.565</u>
d2 or hk0 (Camera K/zero row dist.)	<u>9.13</u>	9.055	<u>8.602</u>	<u>9.508</u>
d1 or hk1 (Camera K/slant vector dist.)	<u>4.56</u>	4.482	<u>4.258</u>	<u>4.706</u>
Ratio of hk0/hk1	<u>2.00</u>	2.02	<u>1.919</u>	<u>2.121</u>
Angle of Slant Vector (Measured)	60.0	60.3	<u>57.314</u>	<u>63.347</u>

From SAED Reference Book, "unknown" diffraction pattern was found to

be that of: Actinolite (Reference JCPDS)

By: EJWP

5209 10/21/2013

With a Zone Axis of: [1 0 0]

Preliminary Identification was: ☒ CORRECT

☐ INCORRECT

percent accuracy to date: _____



AMPHIBOLE SAED INDEXING FORM

EMSL Order ID : Actinolite DATE: 09/27/13

Indexing of negative number: 29035 SCOPE # 27-2

Reference / Sample Number: Actinolite

Preliminary ID: Actinolite By: EJWP

Using Camera Constant of: 22.82 mm Angstroms

Determined from negative number: 29010

Measured Inter-row spacing: 4.37 mm
Mean Distance between spots on Center row (d2): 4.55 mm
Mean distance between spots on slant vector (d1): 8.92 mm

	Calculated	Ref	- 5%	+ 5%
Inter-row spacing (Ångstroms)	<u>5.22</u>	<u>5.300</u>	<u>5.035</u>	<u>5.565</u>
d2 or hk0 (Camera K/zero row dist.)	<u>5.02</u>	5.099	<u>4.844</u>	<u>5.354</u>
d1 or hk1 (Camera K/slant vector dist.)	<u>2.56</u>	2.560	<u>2.432</u>	<u>2.688</u>
Ratio of hk0/hk1	<u>1.96</u>	1.03	<u>0.979</u>	<u>1.082</u>
Angle of Slant Vector (Measured)	<u>85.0</u>	83.0	<u>78.850</u>	<u>87.150</u>

From SAED Reference Book, "unknown" diffraction pattern was found to
be that of: Actinolite (Reference JCPDS) By: EJWP
With a Zone Axis of: [3 -1 -10]

Preliminary Identification was: ☒ CORRECT

☐ INCORRECT

percent accuracy to date: _____



AMPHIBOLE SAED INDEXING FORM

EMSL Order ID : Diopside

DATE: 09/27/13

Indexing of negative number: 29041

SCOPE # 27-2

Reference / Sample Number: Diopside

Preliminary ID: Diopside

By: RSP

Using Camera Constant of: 22.82 mm Angstroms

Determined from negative number: 29010

Measured Inter-row spacing: 4.4 mm
Mean Distance between spots on Center row (d2): 3.63 mm
Mean distance between spots on slant vector (d1): 4.57 mm

	Calculated	Ref	- 5%	+ 5%
Inter-row spacing (Ångstroms)	5.19	5.250	4.988	5.513
d2 or hk0 (Camera K/zero row dist.)	6.29	6.401	6.081	6.721
d1 or hk1 (Camera K/slant vector dist.)	4.99	5.042	4.790	5.294
Ratio of hk0/hk1	1.26	1.27	1.207	1.334
Angle of Slant Vector (Measured)	79.0	78.4	74.480	82.320

From SAED Reference Book, "unknown" diffraction pattern was found to
be that of: Diopside By: EJWP *ESJP*
With a Zone Axis of: [-1 -1 0]

Preliminary Identification was: ☒ CORRECT

☐ INCORRECT

percent accuracy to date: _____



AMPHIBOLE SAED INDEXING FORM

EMSL Order ID : Diopside

DATE: 09/27/13

Indexing of negative number: 29047

SCOPE # 27-2

Reference / Sample Number: DIOPSIDE

Preliminary ID: Diopside

By: EJWP

Using Camera Constant of: 22.82 mm Angstroms

Determined from negative number: 29010

Measured Inter-row spacing: 4.33 mm
Mean Distance between spots on Center row (d2): 6.2 mm
Mean distance between spots on slant vector (d1): 9 mm

	Calculated	Ref	- 5%	+ 5%
Inter-row spacing (Ångstroms)	<u>5.27</u>	<u>5.250</u>	<u>4.988</u>	<u>5.513</u>
d2 or hk0 (Camera K/zero row dist.)	<u>3.68</u>	3.621	<u>3.440</u>	<u>3.802</u>
d1 or hk1 (Camera K/slant vector dist.)	<u>2.54</u>	2.552	<u>2.424</u>	<u>2.680</u>
Ratio of hk0/hk1	<u>1.45</u>	1.42	<u>1.349</u>	<u>1.491</u>
Angle of Slant Vector (Measured)	85.0	84.4	<u>80.180</u>	<u>88.620</u>

From SAED Reference Book, "unknown" diffraction pattern was found to be that of: Diopside

By: EJWP ISUP

With a Zone Axis of: [-4 2 2]

Preliminary Identification was: ☒ X CORRECT

☐ INCORRECT

percent accuracy to date: _____



AMPHIBOLE SAED INDEXING FORM

EMSL Order ID : Diopside

DATE: 09/27/13

Indexing of negative number: 29057

SCOPE # 27-2

Reference / Sample Number: Diopside

Preliminary ID: Diopside

By: EJWP

Using Camera Constant of: 22.82 mm Angstroms

Determined from negative number: 29010

Measured Inter-row spacing:

4.37 mm

Mean Distance between spots on Center row (d2):

3.64 mm

Mean distance between spots on slant vector (d1):

8.92 mm

	Calculated	Ref	- 5%	+ 5%
Inter-row spacing (Ångstroms)	<u>5.22</u>	<u>5.250</u>	<u>4.988</u>	<u>5.513</u>
d2 or hk0 (Camera K/zero row dist.)	<u>6.27</u>	6.401	<u>6.081</u>	<u>6.721</u>
d1 or hk1 (Camera K/slant vector dist.)	<u>2.56</u>	2.550	<u>2.423</u>	<u>2.678</u>
Ratio of hk0/hk1	<u>2.45</u>	2.51	<u>2.385</u>	<u>2.636</u>
Angle of Slant Vector (Measured)	80.0	79.9	<u>75.905</u>	<u>83.895</u>

From SAED Reference Book, "unknown" diffraction pattern was found to be that of: Diopside

By: EJWP *EJWP*

With a Zone Axis of: [2 2 2]

Preliminary Identification was: ☒ CORRECT

☐ INCORRECT

percent accuracy to date: _____



AMPHIBOLE SAED INDEXING FORM

EMSL Order ID : Actinolite

DATE: 9/27/13

Indexing of negative number: 29019

SCOPE # 27-2

Reference / Sample Number: _____

Preliminary ID: Actinolite

By: EJWP

Using Camera Constant of: 22.82 mm Angstroms

Determined from negative number: 29010

Measured Inter-row spacing: 4.4 mm
Mean Distance between spots on Center row (d2): 4.52 mm
Mean distance between spots on slant vector (d1): 4.75 mm

	Calculated	Ref	- 5%	+ 5%
Inter-row spacing (Ångstroms)	<u>5.2</u>	<u>5.3</u>	<u>5.035</u>	<u>5.565</u>
d2 or hk0 (Camera K/zero row dist.)	<u>5.05</u>	5.099	<u>4.844</u>	<u>5.354</u>
d1 or hk1 (Camera K/slant vector dist.)	<u>4.80</u>	4.931	<u>4.684</u>	<u>5.178</u>
Ratio of hk0/hk1	<u>1.05</u>	1.03	<u>0.982</u>	<u>1.086</u>
Angle of Slant Vector (Measured)	<u>69.0</u>	68.0	<u>64.591</u>	<u>71.390</u>

From SAED Reference Book, "unknown" diffraction pattern was found to be that of: Actinolite (JCPDS)

By: K. Colson

With a Zone Axis of: [3 1 2]

Preliminary Identification was: ☒ CORRECT

☐ INCORRECT

percent accuracy to date: _____



AMPHIBOLE SAED INDEXING FORM

EMSL Order ID : Actinolite

DATE: 09/27/13

Indexing of negative number: 29031

SCOPE # 27-2

Reference / Sample Number: Actinolite

Preliminary ID: Actinolite

By: EJWP

Using Camera Constant of: 22.82 mm Angstroms

Determined from negative number: 29010

Measured Inter-row spacing: 4.29 mm
Mean Distance between spots on Center row (d2): 2.5 mm
Mean distance between spots on slant vector (d1): 4.99 mm

	Calculated	Ref	- 5%	+ 5%
Inter-row spacing (Ångstroms)	<u>5.32</u>	<u>5.300</u>	<u>5.035</u>	<u>5.565</u>
d2 or hk0 (Camera K/zero row dist.)	<u>9.13</u>	<u>9.055</u>	<u>8.602</u>	<u>9.508</u>
d1 or hk1 (Camera K/slant vector dist.)	<u>4.57</u>	<u>4.482</u>	<u>4.258</u>	<u>4.706</u>
Ratio of hk0/hk1	<u>2.00</u>	<u>2.02</u>	<u>1.919</u>	<u>2.121</u>
Angle of Slant Vector (Measured)	<u>60.5</u>	<u>60.3</u>	<u>57.314</u>	<u>63.347</u>

From SAED Reference Book, "unknown" diffraction pattern was found to

be that of: Actinolite (JCPDS)

By: K. Colberg

With a Zone Axis of: [1 0 0]

Preliminary Identification was: ☒ CORRECT

☐ INCORRECT

percent accuracy to date: _____



AMPHIBOLE SAED INDEXING FORM

EMSL Order ID : Actinolite

DATE: 9/27/13

Indexing of negative number: 29035

SCOPE # 27-2

Reference / Sample Number: _____

Preliminary ID: Actinolite

By: EJWP

Using Camera Constant of: 22.82 mm Angstroms

Determined from negative number: 29010

Measured Inter-row spacing:

4.38 mm

Mean Distance between spots on Center row (d2):

4.54 mm

Mean distance between spots on slant vector (d1):

8.88 mm

	Calculated	Ref	- 5%	+ 5%
Inter-row spacing (Ångstroms)	<u>5.21</u>	<u>5.300</u>	<u>5.035</u>	<u>5.565</u>
d2 or hk0 (Camera K/zero row dist.)	<u>5.03</u>	<u>5.099</u>	<u>4.844</u>	<u>5.354</u>
d1 or hk1 (Camera K/slant vector dist.)	<u>2.57</u>	<u>2.569</u>	<u>2.441</u>	<u>2.697</u>
Ratio of hk0/hk1	<u>1.96</u>	<u>1.99</u>	<u>1.886</u>	<u>2.084</u>
Angle of Slant Vector (Measured)	<u>83.0</u>	<u>83.0</u>	<u>78.879</u>	<u>87.182</u>

From SAED Reference Book, ^{KC10/22/13} "unknown" diffraction pattern was found to be that of: Actinolite (JCPDS)


By: K Colberg

With a Zone Axis of: [3 -1 -10]

Preliminary Identification was: ☒ CORRECT

☐ INCORRECT

percent accuracy to date: _____

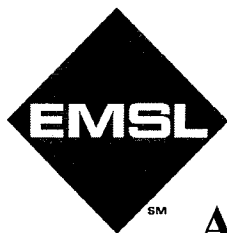

AMPHIBOLE SAED INDEXING FORM

EMSL Order ID : Diopside

DATE: 9/27/13

Indexing of negative number: 29041

SCOPE # 27-2

Reference / Sample Number: _____

Preliminary ID: Diopside

By: RSP

Using Camera Constant of: 22.82 mm Angstroms

Determined from negative number: 29010

Measured Inter-row spacing: 4.45 mm
Mean Distance between spots on Center row (d2): 3.59 mm
Mean distance between spots on slant vector (d1): 4.55 mm

	Calculated	Ref	- 5%	+ 5%
Inter-row spacing (Ångstroms)	<u>5.12</u>	5.250	<u>4.988</u>	<u>5.513</u>
d2 or hk0 (Camera K/zero row dist.)	<u>6.36</u>	6.401	<u>6.081</u>	<u>6.721</u>
d1 or hk1 (Camera K/slant vector dist.)	<u>5.02</u>	5.042	<u>4.790</u>	<u>5.294</u>
Ratio of hk0/hk1	<u>1.27</u>	1.27	<u>1.207</u>	<u>1.334</u>
Angle of Slant Vector (Measured)	<u>79.0</u>	78.4	<u>74.442</u>	<u>82.278</u>

From SAED Reference Book, "unknown" diffraction pattern was found to be that of: Diopside

By: K Colberg

With a Zone Axis of: [-1 -1 0]

Preliminary Identification was: ☒ CORRECT

☐ INCORRECT

percent accuracy to date: _____



AMPHIBOLE SAED INDEXING FORM

EMSL Order ID : Diopside

DATE: 09/27/13

Indexing of negative number: 29047

SCOPE # 27-2

Reference / Sample Number: _____

Preliminary ID: Diopside

By: EJWP

Using Camera Constant of: 22.82 mm Angstroms

Determined from negative number: 29010

Measured Inter-row spacing: 4.4 mm
Mean Distance between spots on Center row (d2): 6.2 mm
Mean distance between spots on slant vector (d1): 8.9 mm

	Calculated	Ref	- 5%	+ 5%
Inter-row spacing (Ångstroms)	5.19	5.250	4.988	5.513
d2 or hk0 (Camera K/zero row dist.)	3.68	3.621	3.440	3.802
d1 or hk1 (Camera K/slant vector dist.)	2.56	2.552	2.424	2.680
Ratio of hk0/hk1	1.44	1.42	1.349	1.491
Angle of Slant Vector (Measured)	84.5	84.4	80.180	88.620

From SAED Reference Book, "unknown" diffraction pattern was found to
be that of: Diopside

By: K Colby

With a Zone Axis of: [-4 2 2]

Preliminary Identification was: ☒ CORRECT

☐ INCORRECT

percent accuracy to date: _____



AMPHIBOLE SAED INDEXING FORM

EMSL Order ID : Diopside

DATE: 9/27/13

Indexing of negative number: 29057

SCOPE # 27-2

Reference / Sample Number: _____

Preliminary ID: Diopside

By: EJWP

Using Camera Constant of: 22.82 mm Angstroms

Determined from negative number: 29010

Measured Inter-row spacing: 4.38 mm
Mean Distance between spots on Center row (d2): 3.51 mm
Mean distance between spots on slant vector (d1): 8.93 mm

	Calculated	Ref	- 5%	+ 5%
Inter-row spacing (Ångstroms)	<u>5.21</u>	<u>5.250</u>	<u>4.988</u>	<u>5.513</u>
d2 or hk0 (Camera K/zero row dist.)	<u>6.50</u>	<u>6.401</u>	<u>6.081</u>	<u>6.721</u>
d1 or hk1 (Camera K/slant vector dist.)	<u>2.56</u>	<u>2.526</u>	<u>2.400</u>	<u>2.652</u>
Ratio of hk0/hk1	<u>2.54</u>	<u>2.53</u>	<u>2.407</u>	<u>2.661</u>
Angle of Slant Vector (Measured)	<u>77.5</u>	<u>77.2</u>	<u>73.378</u>	<u>81.102</u>

From SAED Reference Book, "unknown" diffraction pattern was found to be that of: Diopside

By: K Colberg

With a Zone Axis of: [2 -2 2]

Preliminary Identification was: ☒ CORRECT

☐ INCORRECT

percent accuracy to date: _____



AMPHIBOLE SAED INDEXING FORM

EMSL Order ID : Actinolite

DATE: 09/27/13

Indexing of negative number: 29019

SCOPE # 27-2

Reference / Sample Number: Actinolite

Preliminary ID: Actinolite

By: EJWP

Using Camera Constant of: 22.82 mm Angstroms

Determined from negative number: 29010

Measured Inter-row spacing: 4.5 mm
Mean Distance between spots on Center row (d2): 4.45 mm
Mean distance between spots on slant vector (d1): 4.81 mm

	Calculated	Ref	- 5%	+ 5%
Inter-row spacing (Ångstroms)	<u>5.07</u>	<u>5.30</u>	<u>5.035</u>	<u>5.565</u>
d2 or hk0 (Camera K/zero row dist.)	<u>5.13</u>	5.099	<u>4.844</u>	<u>5.354</u>
d1 or hk1 (Camera K/slant vector dist.)	<u>4.74</u>	4.931	<u>4.684</u>	<u>5.178</u>
Ratio of hk0/hk1	<u>1.08</u>	1.03	<u>0.979</u>	<u>1.082</u>
Angle of Slant Vector (Measured)	68.0	68.0	<u>64.600</u>	<u>71.400</u>

From SAED Reference Book, "unknown" diffraction pattern was found to

be that of: Actinolite

By: RSP

With a Zone Axis of: [3 1 2]

Preliminary Identification was: ☒ CORRECT

☐ INCORRECT

percent accuracy to date: _____



AMPHIBOLE SAED INDEXING FORM

EMSL Order ID : Actinolite

DATE: 09/27/13

Indexing of negative number: 29031

SCOPE # 27-2

Reference / Sample Number: Actinolite

Preliminary ID: Actinolite

By: EJWP

Using Camera Constant of: 22.82 mm Angstroms

Determined from negative number: 29010

Measured Inter-row spacing: 4.25 mm
Mean Distance between spots on Center row (d2): 2.5 mm
Mean distance between spots on slant vector (d1): 5 mm

	Calculated	Ref	- 5%	+ 5%
Inter-row spacing (Ångstroms)	<u>5.37</u>	<u>5.30</u>	<u>5.035</u>	<u>5.565</u>
d2 or hk0 (Camera K/zero row dist.)	<u>9.13</u>	9.055	<u>8.602</u>	<u>9.508</u>
d1 or hk1 (Camera K/slant vector dist.)	<u>4.56</u>	4.482	<u>4.258</u>	<u>4.706</u>
Ratio of hk0/hk1	<u>2.00</u>	2.02	<u>1.919</u>	<u>2.121</u>
Angle of Slant Vector (Measured)	<u>60.0</u>	60.3	<u>57.314</u>	<u>63.347</u>

From SAED Reference Book, "unknown" diffraction pattern was found to

be that of: Actinolite (JCPDS)

By: RSP

With a Zone Axis of: [1 0 0]

Preliminary Identification was: ☒ CORRECT

☐ INCORRECT

percent accuracy to date: _____



AMPHIBOLE SAED INDEXING FORM

EMSL Order ID : Actinolite

DATE: 09/27/13

Indexing of negative number: 29035

SCOPE # 27-2

Reference / Sample Number: Actinolite

Preliminary ID: Actinolite

By: EJWP

Using Camera Constant of: 22.82 mm Angstroms

Determined from negative number: 29010

Measured Inter-row spacing: 4.5 mm
Mean Distance between spots on Center row (d2): 4.56 mm
Mean distance between spots on slant vector (d1): 8.88 mm

	Calculated	Ref	- 5%	+ 5%
Inter-row spacing (Ångstroms)	<u>5.07</u>	<u>5.30</u>	<u>5.035</u>	<u>5.565</u>
d2 or hk0 (Camera K/zero row dist.)	<u>5.00</u>	5.099	<u>4.844</u>	<u>5.354</u>
d1 or hk1 (Camera K/slant vector dist.)	<u>2.57</u>	2.569	<u>2.441</u>	<u>2.697</u>
Ratio of hk0/hk1	<u>1.95</u>	1.99	<u>1.886</u>	<u>2.084</u>
Angle of Slant Vector (Measured)	<u>85.0</u>	83.0	<u>78.850</u>	<u>87.150</u>

From SAED Reference Book, "unknown" diffraction pattern was found to

be that of: Actinolite (JCPDS)

By: RSP

With a Zone Axis of: [3 -1 -10]

Preliminary Identification was: ☒ CORRECT

☐ INCORRECT

percent accuracy to date: _____



AMPHIBOLE SAED INDEXING FORM

EMSL Order ID : Diopside

DATE: 09/27/13

Indexing of negative number: 29041

SCOPE # 27-2

Reference / Sample Number: Diopside

Preliminary ID: Diopside

By: RSP

Using Camera Constant of: 22.82 mm Angstroms

Determined from negative number: 29010

Measured Inter-row spacing: 4.5 mm
Mean Distance between spots on Center row (d2): 3.68 mm
Mean distance between spots on slant vector (d1): 4.63 mm

	Calculated	Ref	- 5%	+ 5%
Inter-row spacing (Ångstroms)	<u>5.07</u>	<u>5.25</u>	<u>4.988</u>	<u>5.513</u>
d2 or hk0 (Camera K/zero row dist.)	<u>6.20</u>	6.401	<u>6.081</u>	<u>6.721</u>
d1 or hk1 (Camera K/slant vector dist.)	<u>4.93</u>	5.042	<u>4.790</u>	<u>5.294</u>
Ratio of hk0/hk1	<u>1.26</u>	1.27	<u>1.207</u>	<u>1.334</u>
Angle of Slant Vector (Measured)	<u>79.0</u>	78.4	<u>74.480</u>	<u>82.320</u>

From SAED Reference Book, "unknown" diffraction pattern was found to

be that of: Diopside

By: RSP

With a Zone Axis of: [-1 -1 0]

Preliminary Identification was: ☒ CORRECT

☐ INCORRECT

percent accuracy to date: _____



AMPHIBOLE SAED INDEXING FORM

EMSL Order ID : Diopside

DATE: 09/27/13

Indexing of negative number: 29047

SCOPE # 27-2

Reference / Sample Number: Diopside

Preliminary ID: Diopside

By: RSP

Using Camera Constant of: 22.82 mm Angstroms

Determined from negative number: 29010

Measured Inter-row spacing: 4.5 mm
Mean Distance between spots on Center row (d2): 6.15 mm
Mean distance between spots on slant vector (d1): 8.87 mm

	Calculated	Ref	- 5%	+ 5%
Inter-row spacing (Ångstroms)	<u>5.07</u>	<u>5.25</u>	<u>4.988</u>	<u>5.513</u>
d2 or hk0 (Camera K/zero row dist.)	<u>3.71</u>	3.621	<u>3.440</u>	<u>3.802</u>
d1 or hk1 (Camera K/slant vector dist.)	<u>2.57</u>	2.552	<u>2.424</u>	<u>2.680</u>
Ratio of hk0/hk1	<u>1.44</u>	1.42	<u>1.349</u>	<u>1.491</u>
Angle of Slant Vector (Measured)	<u>85.0</u>	84.4	<u>80.180</u>	<u>88.620</u>

From SAED Reference Book, "unknown" diffraction pattern was found to be that of: Diopside

By: RSP

With a Zone Axis of: [-4 2 2]

Preliminary Identification was: ☒ CORRECT

☐ INCORRECT

percent accuracy to date: _____



AMPHIBOLE SAED INDEXING FORM

EMSL Order ID : Diopside

DATE: 09/27/13

Indexing of negative number: 29057

SCOPE # 27-2

Reference / Sample Number: Diopside

Preliminary ID: Diopside

By: EJWP

Using Camera Constant of: 22.82 mm Angstroms

Determined from negative number: 29010

Measured Inter-row spacing: 4.5 mm
Mean Distance between spots on Center row (d2): 3.54 mm
Mean distance between spots on slant vector (d1): 9 mm

	Calculated	Ref	- 5%	+ 5%
Inter-row spacing (Ångstroms)	<u>5.07</u>	<u>5.25</u>	<u>4.988</u>	<u>5.512</u>
d2 or hk0 (Camera K/zero row dist.)	<u>6.45</u>	6.401	<u>6.081</u>	<u>6.721</u>
d1 or hk1 (Camera K/slant vector dist.)	<u>2.54</u>	2.550	<u>2.423</u>	<u>2.678</u>
Ratio of hk0/hk1	<u>2.54</u>	2.51	<u>2.385</u>	<u>2.636</u>
Angle of Slant Vector (Measured)	80.0	79.9	<u>75.905</u>	<u>83.895</u>

From SAED Reference Book, "unknown" diffraction pattern was found to

be that of: Diopside

By: RSP

With a Zone Axis of: [2 2 2]

Preliminary Identification was: ☒ CORRECT

☐ INCORRECT

percent accuracy to date: _____

Non-conformance/Corrective Action Report

Section 1: Complete at time problem is identified				Lab: EMSL27	
Department:	Asbestos	CAR# (Assigned by Lab Mgr):	1310-01		
Date Initiated:	10/10/13	Person Initiating CAR:	R. Pescador		
Brief Problem Statement (30 Characters Max)	Documentation Issues				
Full Description of Non-conformance/ problem: (Note things like analyst names, order#, sample#, who reported, etc.)	Internal ILR 2713 271300800-0001. Grid opening and grid information was incorrectly documented.				
Section 2: Complete when initial evaluation is performed					
Evaluation/ Investigation Details	From R. Denton: For Sample 271300800-0002 1 LA structure was listed on the benchsheet and Verified GO Map in GO C2: C8. The EDXA spectra you provided said that the structure was located in GO C1:G8. This structure was not located during the interlab verified analysis in either grid opening listed on the bench sheet/verified map or EDXA. During the reconciliation, the structure was located in GO C2:G8. This improper documentation lead to the false positive. Improper documentation. Labelling spectra based on the information collected from the benchsheet. The grid being analyzed (C2) was not recorded on the benchsheet at the time the spectrum as named and saved. I documented the grid on the spectrum to be C1 instead of C2 since C2 was not on the benchsheet at that time. GO C8 was recorded instead G8 du to the grid being oriented in a way that "G" looked like a "C".				
Evaluating Party & Date Completed:	R. Denton 10/10/13				
Remedial Actions Taken:	Have another personnel review paperwork before it is sent out or reported.				
Is Root Cause & Corrective Action Required (Yes/No)?	Yes	If not, provide justification for this decision			
Client notification required? <i>If yes, record contact details (who, when, how)?</i>	No				
Is Work Stop Necessary? <i>If yes, record details (who, when, until when)?</i>	No				
If Work Stop, record details of when it was lifted and who authorized:					
Section 3: Complete whenever a corrective action is required					
Root Cause:	Improper documentation was caused by not completing analytical worksheet correctly prior to moving to the next grid opening.				
Corrective Action:	When saving spectra and documenting GO's, switch to low mag and verify the correct opening is being documented. Analyst will ensure that directly upon moving to a new Grid, the grid ID will be documented before proceeding to a grid opening. Prior to recording the Grid opening the analyst will switch to low mag and ensure the GO ID. Once GO ID is confirmed the analyst will document it on the analytical worksheet.				
Corrective Action due date:	10/10/13	Assigned to:	R. Pescador		
Evidence of commitment and compliance: (e.g., revised documents, new documents, corrected records, etc. Docs can be inserted using Insert>Object)	An extra level of QA will be performed at the lab level. There shall be 2 initials and date on the first page of the ICOC adjacent to the due date/time to show that at least 2 personnel have reviewed the paperwork.				
Manager Signoff:	R. Pescador	Date:	10/10/13		
Section 4: Complete after implementation of corrective action					
Follow-Up Due Date:	10/10/13	Follow-up Notes (i.e., how was effectiveness determined, findings):			
Follow-Up Completed By:	E. Wyatt-Pescador				
Date Completed:	10/10/13	Outcome:	Effective/ No Recurrences	New CAR # (when necessary):	

Matched Structures per GO

EMSL Order ID		271300573		Sample ID:		SM-10080		Original Lab/QC Lab		EMSL27/EMSL04		Printed: 10/29/2013		
Grid Opening ID	Re-preparation					Interlab					NISTIR 5351 QC CODE			Comments
Grid:GO	Total Str #	Str Type	Length	Width	Mineral ID	Total Str #	Str Type	Length	Width	Mineral ID	1st Read	2nd Read	RECONCILE	
1	E1 : B2	ND					ND							
2	E1 : B4	ND					ND							
3	E1 : B6	ND				0	F	2.6	0.5	LA	TN	TN		SAP references modification that could not be determined or found on eRoom. Looking through other available documents, it was concluded that analysis was PCME >5.0um.
4	E1 : B8	ND					ND							
5	E1 : B10	ND					ND							
6	E2 : E1	ND					ND							
7	E2 : E3	ND					ND							
8	E2 : E5	ND					ND							
9	E2 : E7	ND					ND							
10	E2 : E9	ND					ND							
11	E3 : J1	ND					ND							
12	E3 : J3	ND					ND							
13	E3 : J5	ND					ND							
14	E3 : J7	ND					ND							
15	E3 : J9	ND					ND							

Matched Structures per GO for Reconciliation Analysis													Printed: 10/29/2013	
EMSL Order ID		271300675		Sample ID:		P3-00774		Original Lab/QC Lab		EMSL27/EMSL04				
Grid Opening ID	Re-preparation					Interlab					Results		Comments	
Grid:GO	Total Str #	Str Type	Length	Width	Mineral ID	Total Str #	Str Type	Length	Width	Mineral ID	Conclusion	Failures due to		
1	G1 : A1		ND				ND				Pass			
2	G1 : A2		ND				ND				Pass			
3	G1 : A3		ND				ND				Pass			
4	G1 : A4		ND				ND				Pass			
5	G1 : A5		ND				ND				Pass			
6	G1 : A6		ND				ND				Pass			
7	G1 : A7		ND				ND				Pass			
8	G1 : A8		ND				ND				Pass			
9	G1 : A9		ND				ND				Pass			
10	G1 : A10		ND				ND				Pass			
11	G1 : C1		ND				ND				Pass			
12	G1 : C2		ND				ND				Pass			
13	G1 : C3		ND				ND				Pass			
14	G1 : C4		ND				ND				Pass			

Reconciliation Analyst: _____

Date: _____

Scope: _____

Comments:

Matched Structures per GO for Reconciliation Analysis													Printed: 10/29/2013	
EMSL Order ID		271300675		Sample ID:		P3-00774		Original Lab/QC Lab		EMSL27/EMSL04				
Grid Opening ID	Re-preparation					Interlab					Results		Comments	
Grid:GO	Total Str #	Str Type	Length	Width	Mineral ID	Total Str #	Str Type	Length	Width	Mineral ID	Conclusion	Failures due to		
15	G1 : C5	ND					ND				Pass			
16	G1 : C6	ND					ND				Pass		GRID TORN	
17	G1 : C7	ND					ND				Pass		GRID TORN	
18	G1 : C8	ND					ND				Pass			
19	G1 : C9	ND					ND				Pass			
20	G1 : C10	ND					ND				Pass			
21	G1 : E1	ND					ND				Pass			
22	G1 : E2	ND					ND				Pass			
23	G1 : E3	ND					ND				Pass			
24	G1 : E4	ND					ND				Pass			
25	G1 : E5	ND					ND				Pass			
26	G1 : E6	ND					ND				Pass			
27	G1 : E7	ND					ND				Pass			
28	G1 : E8	ND					ND				Pass			

Reconciliation Analyst: _____

Date: _____

Scope: _____

Matched Structures per GO for Reconciliation Analysis													Printed: 10/29/2013	
EMSL Order ID		271300675		Sample ID:		P3-00774		Original Lab/QC Lab		EMSL27/EMSL04				
Grid Opening ID	Re-preparation					Interlab					Results		Comments	
Grid:GO	Total Str #	Str Type	Length	Width	Mineral ID	Total Str #	Str Type	Length	Width	Mineral ID	Conclusion	Failures due to		
29	G1 : E9	ND					ND				Pass			
30	G1 : E10	ND					ND				Pass			
31	G1 : G1	ND					ND				Pass			
32	G1 : G2	ND					ND				Pass			
33	G1 : G3	ND					ND				Pass			
34	G1 : G4	ND					ND				Pass			
35	G1 : G5	ND					ND				Pass			
36	G1 : G6	ND					ND				Pass			
37	G1 : G7	ND					ND				Pass			
38	G1 : G8	ND					ND				Pass			
39	G1 : G9	ND					ND				Pass			
40	G1 : G10	ND					ND				Pass			
41	G1 : I1	ND					ND				Pass			
42	G1 : I2	ND					ND				Pass			

Reconciliation Analyst: _____

Date: _____

Scope: _____

Matched Structures per GO for Reconciliation Analysis													Printed: 10/29/2013	
EMSL Order ID		271300675		Sample ID:		P3-00774		Original Lab/QC Lab		EMSL27/EMSL04				
Grid Opening ID	Re-preparation					Interlab					Results		Comments	
Grid:GO	Total Str #	Str Type	Length	Width	Mineral ID	Total Str #	Str Type	Length	Width	Mineral ID	Conclusion	Failures due to		
3	G1 : I3	ND					ND				Pass			
4	G1 : I4	ND					ND				Pass			
5	G1 : I5	ND					ND				Pass			
6	G1 : I6	ND					ND				Pass			
7	G1 : I7	ND					ND				Pass			
8	G1 : I8	ND					ND				Pass			
9	G1 : I9	ND					ND				Pass			
0	G1 : I10	ND					ND				Pass			
1	G2 : B1	ND					ND				Pass			
2	G2 : B2	ND					ND				Pass			
3	G2 : B3	ND					ND				Pass			
4	G2 : B4	ND					ND				Pass			
5	G2 : B5	ND					ND				Pass			
6	G2 : B6	ND					ND				Pass			

Reconciliation Analyst: _____

Date: _____

Scope: _____

Matched Structures per GO for Reconciliation Analysis													Printed: 10/29/2013	
EMSL Order ID		271300675		Sample ID:		P3-00774		Original Lab/QC Lab		EMSL27/EMSL04				
Grid Opening ID	Re-preparation					Interlab					Results		Comments	
Grid:GO	Total Str #	Str Type	Length	Width	Mineral ID	Total Str #	Str Type	Length	Width	Mineral ID	Conclusion	Failures due to		
7 G2 : B7		ND					ND				Pass			
8 G2 : B8		ND					ND				Pass			
9 G2 : B9		ND					ND				Pass			
0 G2 : B10		ND					ND				Pass			
1 G2 : D1		ND					ND				Pass			
2 G2 : D2		ND					ND				Pass			
3 G2 : D3		ND					ND				Pass			
4 G2 : D4		ND					ND				Pass			
5 G2 : D5		ND					ND				Pass			
5 G2 : D6		ND					ND				Pass			
7 G2 : D7		ND					ND				Pass			
8 G2 : D8		ND					ND				Pass			
9 G2 : D9		ND					ND				Pass			
0 G2 : D10		ND					ND				Pass			

Reconciliation Analyst: _____

Date: _____

Scope: _____

Matched Structures per GO for Reconciliation Analysis													Printed: 10/29/2013	
EMSL Order ID		271300675		Sample ID:		P3-00774		Original Lab/QC Lab		EMSL27/EMSL04				
Grid Opening ID	Re-preparation					Interlab					Results		Comments	
Grid:GO	Total Str #	Str Type	Length	Width	Mineral ID	Total Str #	Str Type	Length	Width	Mineral ID	Conclusion	Failures due to		
71	G2 : F1	ND					ND				Pass			
72	G2 : F2	ND					ND				Pass			
73	G2 : F3	ND					ND				Pass			
74	G2 : F4	ND					ND				Pass			
75	G2 : F5	ND					ND				Pass			
76	G2 : F6	ND					ND				Pass			
77	G2 : F7	ND					ND				Pass			
78	G2 : F8	ND					ND				Pass			
79	G2 : F9	ND					ND				Pass			
80	G2 : F10	ND					ND				Pass			
81	G2 : H1	ND					ND				Pass			
82	G2 : H2	ND					ND				Pass			
83	G2 : H3	ND					ND				Pass			
84	G2 : H4	ND					ND				Pass			

Reconciliation Analyst: _____

Date: _____

Scope: _____

Matched Structures per GO for Reconciliation Analysis													Printed: 10/29/2013		
EMSL Order ID			271300675		Sample ID:		P3-00774		Original Lab/QC Lab		EMSL27/EMSL04				
Grid Opening ID	Re-preparation					Interlab					Results		Comments		
Grid:GO	Total Str #	Str Type	Length	Width	Mineral ID	Total Str #	Str Type	Length	Width	Mineral ID	Conclusion	Failures due to			
5 G2 : H5		ND					ND				Pass				
5 G2 : H6		ND					ND				Pass				
7 G2 : H7		ND					ND				Pass				
8 G2 : H8		ND					ND				Pass				
9 G2 : H9		ND					ND				Pass				
0 G2 : H10		ND					ND				Pass				
1 G2 : J1		ND					ND				Pass				
2 G2 : J2		ND					ND				Pass				
3 G2 : J3		ND					ND				Pass				
4 G2 : J4		ND					ND				Pass				
5 G2 : J5		ND					ND				Pass				
5 G2 : J6		ND					ND				Pass		GRID TORN		
7 G2 : J7		ND					ND				Pass		GRID TORN		
8 G2 : J8		ND					ND				Pass				

Reconciliation Analyst: _____

Date: _____

Scope: _____

Matched Structures per GO for Reconciliation Analysis

Printed: 10/29/2013

Reconciliation Analyst: _____

Date: _____

Scope: _____

EMSL Order ID 271300675 Sample ID: P3-00774 Original Lab/QC Lab EMSL27/EMSL04														
Grid Opening ID	Re-preparation					Interlab					Results		Comments	
Grid:GO	Total Str #	Str Type	Length	Width	Mineral ID	Total Str #	Str Type	Length	Width	Mineral ID	Conclusion	Failures due to		
99	G2 : J9		ND									Pass		
100	G2 : J10		ND									Pass		
101	H1 : B1		ND									Pass		
102	H1 : B2		ND									Pass		
103	H1 : B3		ND									Pass		
104	H1 : B4		ND									Pass		
105	H1 : B5		ND									Pass		
106	H1 : B6		ND									Pass		
107	H1 : B7		ND									Pass		
108	H1 : B8		ND									Pass		
109	H1 : B9		ND									Pass		
110	H1 : B10		ND									Pass		
111	H1 : C1		ND									Pass		
112	H1 : C2		ND									Pass		

Matched Structures per GO for Reconciliation Analysis

Printed: 10/29/2013

Reconciliation Analyst: _____

Date: _____

Scope: _____

Matched Structures per GO for Reconciliation Analysis														Printed: 10/29/2013	
EMSL Order ID 271300675						Sample ID: P3-00774		Original Lab/QC Lab EMSL27/EMSL04							
Grid Opening ID		Re-preparation					Interlab					Results		Comments	
Grid:GO		Total Str #	Str Type	Length	Width	Mineral ID	Total Str #	Str Type	Length	Width	Mineral ID	Conclusion	Failures due to		
113	H1 : C3		ND					ND				Pass			
114	H1 : C4		ND					ND				Pass			
115	H1 : C5		ND					ND				Pass			
116	H1 : C6		ND					ND				Pass			
117	H1 : C7		ND					ND				Pass			
118	H1 : C8		ND					ND				Pass			
119	H1 : C9		ND					ND				Pass			
120	H1 : C10		ND					ND				Pass			
121	H1 : E1		ND					ND				Pass			
122	H1 : E2		ND					ND				Pass			
123	H1 : E3		ND					ND				Pass			
124	H1 : E4		ND					ND				Pass			
125	H1 : E5		ND					ND				Pass			
126	H1 : E6		ND					ND				Pass			

Matched Structures per GO for Reconciliation Analysis

Printed: 10/29/2013

Reconciliation Analyst: _____

Date: _____

Scope: _____

Matched Structures per GO for Reconciliation Analysis													
EMSL Order ID 271300675						Sample ID: P3-00774			Original Lab/QC Lab EMSL27/EMSL04				
Grid Opening ID	Re-preparation					Interlab					Results		Comments
Grid:GO	Total Str #	Str Type	Length	Width	Mineral ID	Total Str #	Str Type	Length	Width	Mineral ID	Conclusion	Failures due to	
127 H1 : E7		ND					ND				Pass		GRID TORN
128 H1 : E8		ND					ND				Pass		
129 H1 : E9		ND					ND				Pass		
130 H1 : E10		ND					ND				Pass		
131 H1 : F1		ND					ND				Pass		
132 H1 : F2		ND					ND				Pass		
133 H1 : F3		ND					ND				Pass		
134 H1 : F4		ND					ND				Pass		
135 H1 : F5		ND					ND				Pass		
136 H1 : F6		ND					ND				Pass		
137 H1 : F7		ND					ND				Pass		
138 H1 : F8		ND					ND				Pass		
139 H1 : F9		ND					ND				Pass		
140 H1 : F10		ND					ND				Pass		

Matched Structures per GO for Reconciliation Analysis

Printed: 10/29/2013

Reconciliation Analyst: _____

Date: _____

Scope: _____

Matched Structures per GO for Reconciliation Analysis													
EMSL Order ID 271300675						Sample ID: P3-00774			Original Lab/QC Lab EMSL27/EMSL04				
Grid Opening ID	Re-preparation					Interlab					Results		Comments
Grid:GO	Total Str #	Str Type	Length	Width	Mineral ID	Total Str #	Str Type	Length	Width	Mineral ID	Conclusion	Failures due to	
141 H1 : G1		ND					ND				Pass		
142 H1 : G2		ND					ND				Pass		
143 H1 : G3		ND					ND				Pass		
144 H1 : G4		ND					ND				Pass		
145 H1 : G5		ND					ND				Pass		
146 H1 : G6		ND					ND				Pass		
147 H1 : G7		ND					ND				Pass		
148 H1 : G8	0	F	4.05	0.3	LA		ND				Fail	Length, Mineral	
149 H1 : G9		ND					ND				Pass		
150 H1 : G10		ND					ND				Pass		
151 H1 : H1		ND					ND				Pass		
152 H1 : H2		ND					ND				Pass		
153 H1 : H3		ND					ND				Pass		
154 H1 : H4		ND					ND				Pass		

Matched Structures per GO for Reconciliation Analysis

Printed: 10/29/2013

Reconciliation Analyst: _____

Date: _____

Scope: _____

Matched Structures per GO for Reconciliation Analysis													
EMSL Order ID 271300675						Sample ID: P3-00774			Original Lab/QC Lab EMSL27/EMSL04				
Grid Opening ID	Re-preparation					Interlab					Results		Comments
Grid:GO	Total Str #	Str Type	Length	Width	Mineral ID	Total Str #	Str Type	Length	Width	Mineral ID	Conclusion	Failures due to	
155 H1 : H5		ND					ND				Pass		
156 H1 : H6		ND					ND				Pass		
157 H1 : H7		ND					ND				Pass		
158 H1 : H8		ND					ND				Pass		
159 H1 : H9		ND					ND				Pass		
160 H1 : H10		ND					ND				Pass		

[illegible]

[illegible]