


INDIAN HEAD DENTAL LAB, INDIAN HEAD, MD

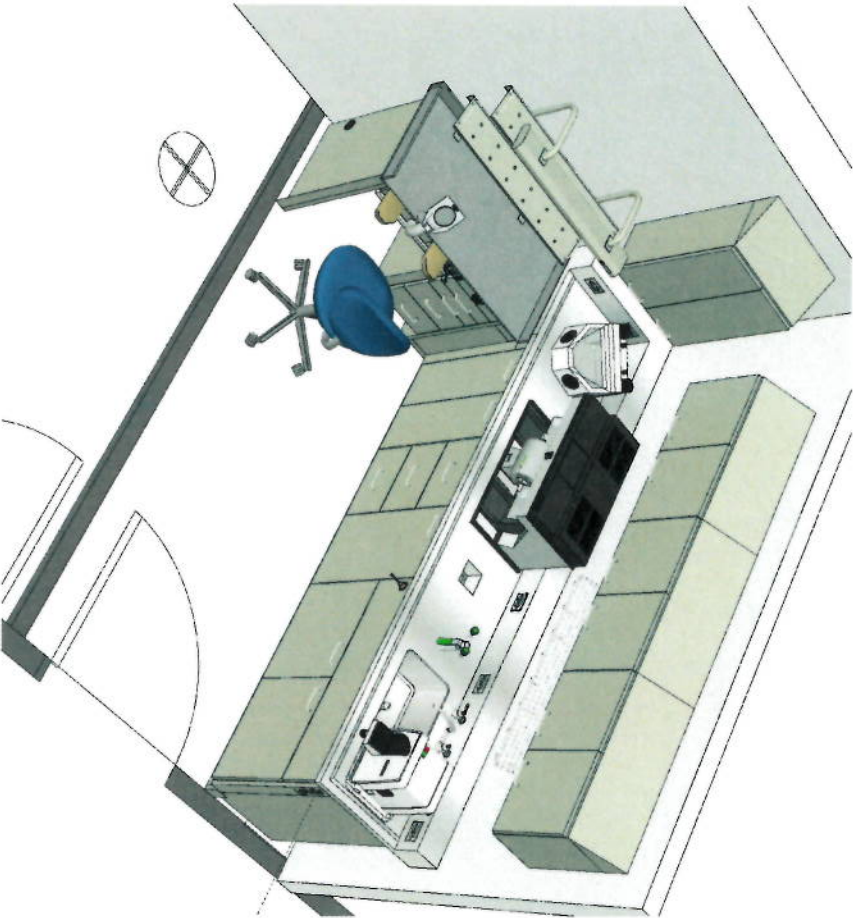

GENERAL NOTES

1. STANDARD BASE CABINETS ARE AT 35 5/8" HEIGHT WITH NO TOP.
2. COUNTER TOP MATERIAL: 1/2" SOLID SURFACE (IN "SANDED VERMILLION" OR "SANDED GRAY"), 1 1/2" STAINLESS STEEL, & 1 1/4" SOAP STONE.
- OTHER COUNTER TOP MATERIAL OPTIONS AVAILABLE UPON REQUEST.
3. ALL WALL HUNG ITEMS EITHER BY NEVINLABS OR BY OTHERS REQUIRE WALL BACKING TO BE INSTALLED BY OTHERS.
4. ALL PIPING, WIRING, CONNECTIONS, SHUT OFF VALVES & FINAL HOOK-UPS TO BE PROVIDED BY OTHERS.
5. NEVINLABS WORK STATION BENCHES, ERGOTECH, MEDALIST, FREEDOM, AND ERGOEZ, ARE PRE-PIPED AND PRE-WIRED, SHIPPED PARTIALLY ASSEMBLED FOR EASY INSTALLATION. SHUT OFF VALVES & FINAL HOOK-UPS TO BE PROVIDED BY OTHERS.
6. PEDESTAL DUST COLLECTOR WORKSTATIONS, IE. WSA40NT, WSA42NT, & WSA4 ARE NOT PRE-PIPED NOR PRE-WIRED. PIPING AND WIRING AND FINAL HOOK-UPS BY OTHERS.
7. EQUIPMENT NOT BY NEVINLABS: CASTING ARM, DRY MODEL TRIMMER, VACUUM, MIXER, COMPUTER, CNC MILLING MACHINE, FURNACE OVEN, ETC. UNLESS NOTED OTHERWISE.
8. ALL UTILITIES IN WALLS, ABOVE AND BELOW NEVIN CASEWORK, ARE BY OTHERS.
9. ELECTRIC OUTLETS IN WALLS ARE NOT FURNISHED BY NEVIN. THEY ARE SHOWN AT A HEIGHT TO CLEAR BACK SPLASH OR BACK SHELF. ARCHITECT OR CONTRACTOR TO COORDINATE QUANTITY AND LOCATION OF ALL ELECTRIC OUTLETS IN WALLS.
10. ALL PIPING, WIRING AND FINAL HOOK UPS ARE THE RESPONSIBILITY OF THE OWNER OR OWNER'S AGENT.
11. ALL DIMENSIONS ARE BASED ON CADFILES OR A HAND-DRAWN SKETCH OR ORIGINAL BLUEPRINT PROVIDED TO NEVINLABS. AS-BUILT DIMENSIONS MUST BE JOB-SITE VERIFIED BY A NEVINLABS DENTAL REPRESENTATIVE.
12. ALL ROOM DIMENSIONS ARE TO BE VERIFIED IN THE FIELD.

THIS DRAWING WAS PREPARED EXCLUSIVELY FOR YOU WITH THE UNDERSTANDING THAT NEVINLABS RESERVES THE SOLE RIGHTS TO ITS POSSESSION AND USE, ANY REPRODUCTION OR USE OF THIS DRAWING THAT IT IS DETRIMENTAL TO OUR INTERESTS IS PROHIBITED.

PROJECT INFORMATION

Naval Branch Health Clinic Building 1600  
Indian Head, MD 20640

JOB CONTACT

CONTACT: Mark Beaupre  
PHONE: 317-413-4556  
EMAIL: mark.beaupre@pattersondental.com  
FAX:

DEALER/CONTRACTOR

COMPANY: PATTERSON DENTAL  
CONTACT: Mark Beaupre  
PHONE: 317-413-4556  
EMAIL: mark.beaupre@pattersondental.com

NEVINLABS CONTACTS

SALES ASSOCIATES

CONTACT: Karl Odden  
PHONE: 484-401-0991  
EMAIL: koddenn@dentaliez.com

FACILITY PLANNER

CONTACT: Anna Jackiewicz  
PHONE: 1-773-624-4330 ext 8211  
EMAIL: ajackiewicz@dentaliez.com

PROJECT MANAGER

CONTACT:  
PHONE:  
EMAIL:

DRAWING INFORMATION

FIXTURE LEGEND

- GAS/AIR/VACUUM - LEVER BALL VALVE, NEEDLE NOSE (206G/206A/206V) 3/8" NPT Final connection by owner
- OXYGEN/GAS/AIR/VACUUM - NEEDLE CONTROL MALE VALVE, NEEDLE NOSE (207S/207G/207A/207V) 3/8" NPT (FINAL CONNECTION BY OWNER)
- AIR QUICK DISCONNECT (209T) 3/8" NPT (FINAL CONNECTION BY OWNER)
- COLD WATER VALVE, 20800 CWV, 3/8" NPT (FINAL CONNECTION BY OWNER)

PROJECT DRAWING INDEX

SHEET NO.	TITLE	SHEET NO.	TITLE
1	PLAN 1/2"=1'-0"	11	
2	ELEVATION 1, TYPICAL DETAILS AND VIEW & UTILITIES	12	
3	ELEVATION 2, TYPICAL DETAILS AND VIEW & UTILITIES	13	
4		14	
5		15	
6		16	
7		17	
8		18	
9		19	
10		20	

REVISIONS

DATE	DESCRIPTION	DATE	DESCRIPTION
A 03/22/2018	dims of the space verified and added utilities	G 00/00/2018	
B		H	
C		J	
D		K	
E		L	
F		M	




5000 SOUTH HALSTED STREET CHICAGO, IL 60609  
PHONE (773) 624-4330  
FAX (773) 624-7337

ISSUE DATE: 01-23-2018  
CAD FILE: 218006  
ORACLE NUMBER: X  
QUOTATION: 20180124012001084339

Attachment A



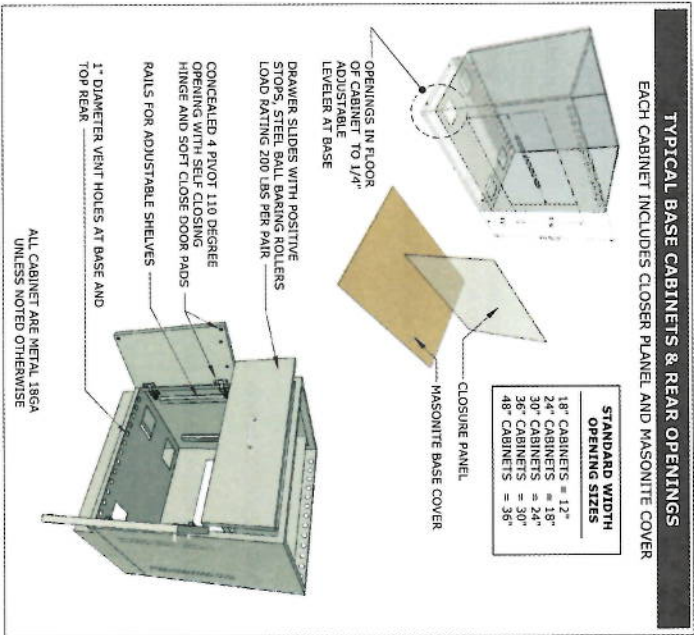

**FINISH MATERIALS COLOR OPTIONS**

**STANDARD METAL FINISHES COLOR OPTIONS**

Blue	Black	White	Red	Green
Yellow	Grey	Light Blue	Dark Blue	Light Green
Light Yellow	Light Grey	Light Blue	Dark Blue	Light Green

**SOLID SURFACE FINISHES COLOR OPTIONS**

White	Black	White	Red	Green
Yellow	Grey	Light Blue	Dark Blue	Light Green
Light Yellow	Light Grey	Light Blue	Dark Blue	Light Green

NOTE: COLORS MAY VARY FROM EXACT FINISHES AND ARE FOR GENERAL REFERENCE ONLY. PLEASE CONTACT YOUR SALES REPRESENTATIVE OR CALL FOR FINISHES FOR PROPOSAL. SAMPLES TO BE SENT TO YOU.

CUSTOM COLORS AVAILABLE UPON REQUEST FOR AN ADDITIONAL CHARGE

**FINISH MATERIAL SELECTION**

**CABINET COLORS SPECIFICATION:**

CABINET BODIES AND FRAMES:	
DOOR AND DRAWER FRONTS:	
OVERHEAD LIGHT HOUSINGS:	
LIGHT SUPPORT:	
LIGHT SHELVES:	

**SOLID SURFACE COUNTER TOPS FINISHES**

SANDED VERMILLION	SANDED GREY
-------------------	-------------

**CUSTOMER PREPARED WORKSTATION HEIGHT**

32"	34"
-----	-----

**STANDARD FINISHES**


**DRAWER PULLS - BRUSHED ALUMINUM FINISHED**

**STANDARD COUNTER TOPS FINISHES**

1 1/2" STAINLESS STEEL - WITH MARINE EDGES

1 1/4" SOAP STONE

PLEASE CONTACT YOUR SALES REPRESENTATIVE OR CALL FOR FINISHES FOR PROPOSAL. SAMPLES TO BE SENT TO YOU.

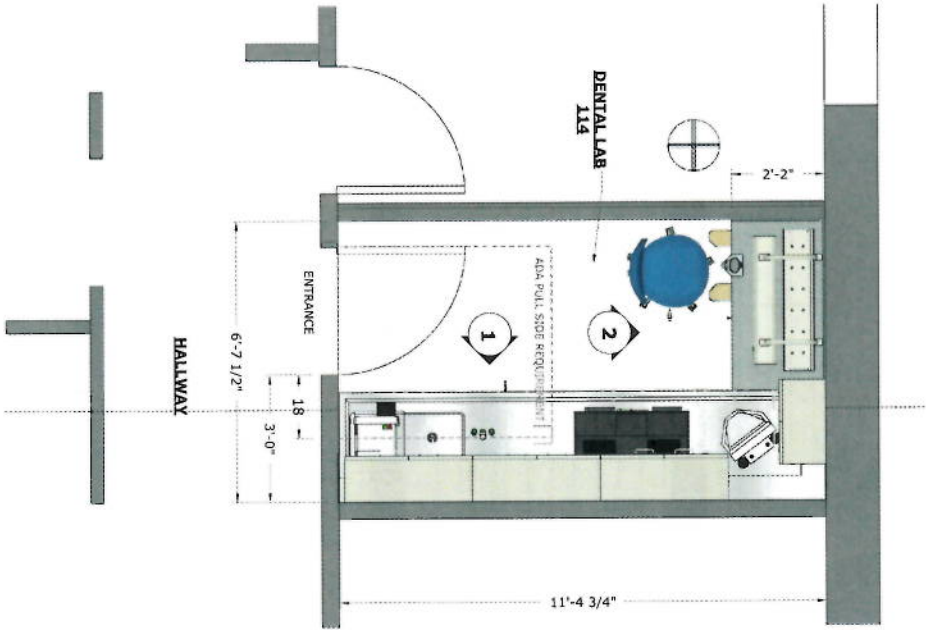

**ROOM DIMENSION AND CEILING HEIGHT TO BE VERIFIED BY OTHERS.**  
ELECTRICAL GAS, AIR, AND UTILITIES ABOVE BACKSLASHES IN WALL AS NEEDED, BY CONTRACTOR

SHEET  
NUMBER

01

LOCATION: Naval Branch Health Clinic Building 1600  
Indian Head, MD 20640

DEALER: PATTERSON DENTAL  
Mark Beaupre  
317-413-4556  
mark.beaupre@pattersondental.com

CAD FILE: 218006

Oracle Number:

DESIGNER: AJ



5000 SOUTH HALSTED STREET CHICAGO, IL 60609  
PHONE (773) 624-4330 FAX (773) 624-7337

INDIAN HEAD DENTAL LAB, INDIAN HEAD, MD


REVISIONS	
DATE	DISCRIPTION
A	00/00/2016
B	
C	
D	
E	

D

NADL  
National Association of  
Dental Laboratories



453DBLSK  
INTEGRALLY WELDED STAINLESS STEEL SINK AND TRAY  
WITH TYPICAL CABINETRY SET UP

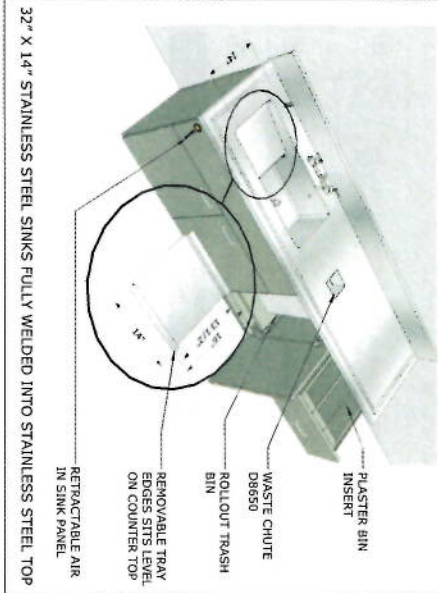


80093  
MODEL TRIMMER WITH BRAKE



- Warning light at base of unit for people with impaired hearing or for noisy work areas. Includes model tray that fits on top of unit.
- 1600 W and 3000 rpm motor.
- Work window: 121 mm (max. of its category)
- Abrasive diameter: 250 mm
- Power: 1600 W
- Speed: 3000 rpm
- Supply: 230 V, 50/60 Hz
- Height: 340 mm
- Width: 420 mm
- Depth: 345 mm
- Weight: 14.5 kg

365 PLASTER BIN  
DRAWER INSERT

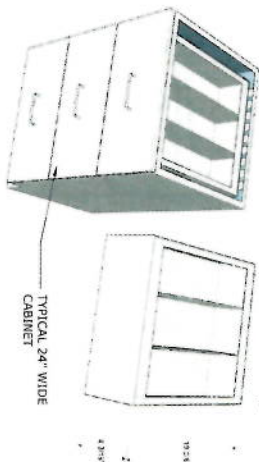


3. COMPARTMENT STORAGE BIN, INSERT
- Overall size 20" wide 19 3/8" long and 8 3/16" deep.
  - Each compartment is 6 11/16" wide by 19 3/8" long.
  - Unit holds up to 150 LBS.

YANIMAN SANDVAC  
DUST COLLECTION SYSTEM



- Voltage 115.00 VAC
- Amperage 5.00 Amps
- Available in 100V, 220V, 230V, 240V @ 50/60Hz


80502 BENCH TOP POLISHING HOOD  
TABLE TOP UNIT FOOTPRINT: 19 3/4" x 31 1/2" x 19 3/4"  
115V - 5.5 AMPS

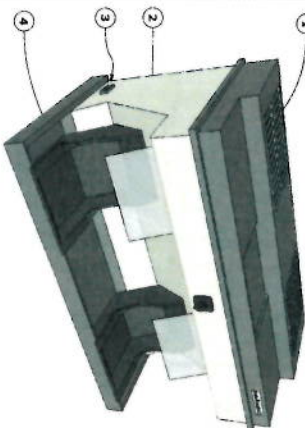

1. (2) PICS NON-WOVEN DISPOSABLE FILTER BAGS UNDER ID
  2. CLEAR SAFETY SHIELDS FOR EACH SPINDLE (ROTATES UP/DOWN)
  3. SAFETY RESET SWITCH (MUST RESET AFTER POWER DISRUPTED)
  4. (2) REMOVABLE RUBBER SPLASH PANS (GAUDES) WITH LIGHT OVER EACH SPINDLE
- CONSTRUCTION  
RUST-PROOF POLYCARBONATE OVER HEAVY GAUGE STEEL FRAME

LATHE SOLD SEPARATELY - NEVINLABS #380WCT LATHE - APPROX. 4AMPS


WIRE SHELVES DETAIL

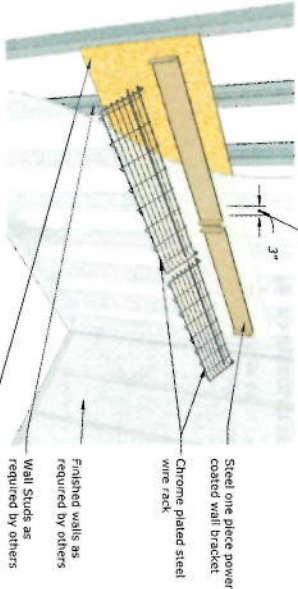

NOTE:  
WIRE SHELVES MOUNTING WALL BRACKETS SIDE BY SIDE, ALLOW MINIMUM 3" BETWEEN FOR ACCESS TO 1/4-20 BOLTS

- Complete 2232250JB unit size: 48" wide x 9" deep x 4" high  
Complete 223223WB unit size: 42" wide x 9" deep x 4" high
- Mounting to stud wall or wall backing: locate center of studs and mark position on walls.
  - Position wall brackets to pick up at least 2 studs and mark the center location.
  - Drill holes in brackets to match the size of the mounting screws you will be using.
  - Mount the brackets to the wall.
  - Check with local hardware store for the correct anchors to use.
- Once the wall bracket is securely fastened to the wall position the wire rack in to the bracket and insert the 1/4-20 bolt through the pre-punched hole from the outside of the bracket. Install the 1/4" nut on the inside of the bracket and tighten.

UTILITY LEGEND


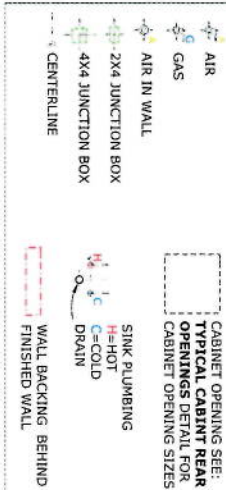


NOTES:

1. SIZE AND QUANTITIES OF JUNCTION BOXES AND OUTLET ARE SUGGESTED.
  2. ACTUAL SIZE AND QUANTITIES NEEDED IN FIELD DETERMINED BY OTHERS.
  3. LOCATIONS OF JUNCTION BOXES AS SHOWN ARE TO CENTERLINE AND TOP OF REAR OPENING OF CABINETS.
  4. UTILITIES IN WALL BY CONTRACTOR (AIR, GAS, ELECTRIC) CONTRACTOR SUPPLY SHUT OFF VALVES AND EXCLUSIVE ONE 120V 20 AMP CIRCUIT PER SINGLE STATION.
- REQUIRED UTILITIES IN WALL BY CONTRACTOR FOR WORKSTATIONS (97418):
- 120V 20 AMP EXCLUSIVE CIRCUIT PER BENCH
  - AIR SHUT OFF WITH 3/8" FEMALE CONNECTION
  - GAS SHUT OFF WITH 1/2" FEMALE CONNECTION
- A = AIR  
G = GAS



ELECTRICAL AND AIR IN WALL ABOVE COUNTER TOP BACK SPLASH AS REQUIRED, BY OTHERS

3D VIEW - POSITION - 1 and 2  
SCALE: N/A

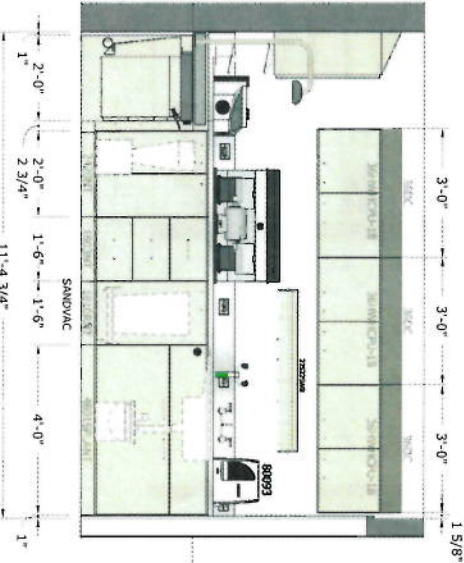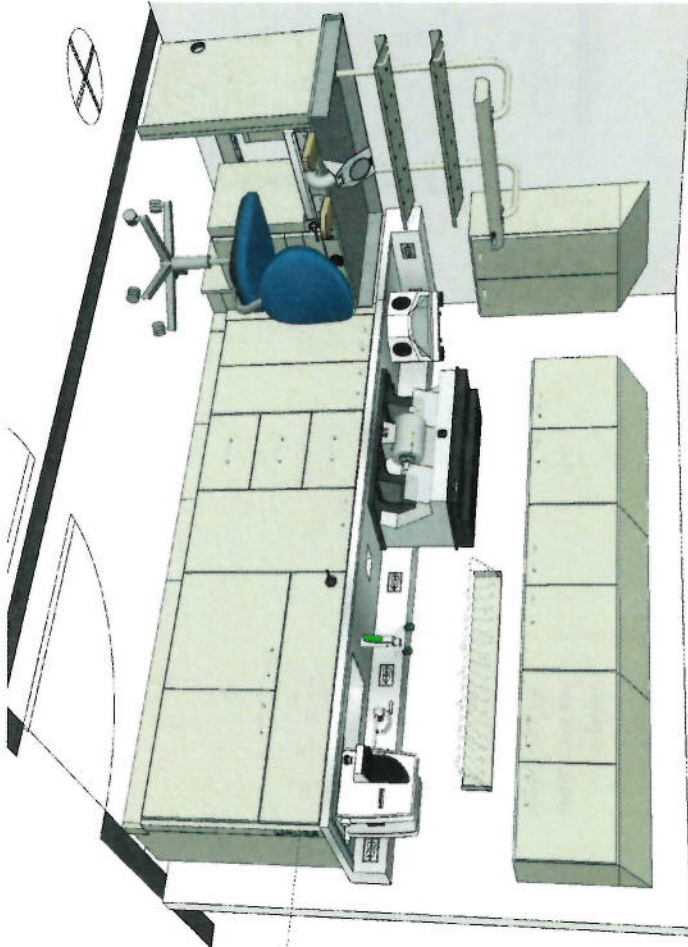

POS 1: SINK AREA  
TOP 1 1/2" STAINLESS STEEL

POS 1: SECTION  
TOP 1 1/2" STAINLESS STEEL

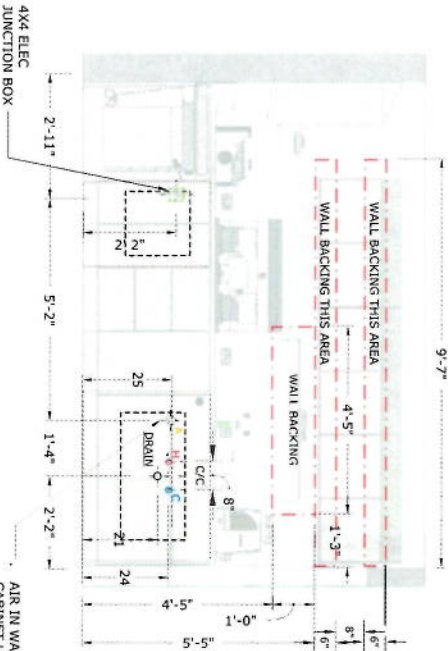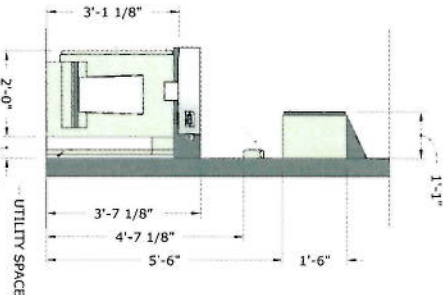

POS 1: SINK AREA  
UTILITIES

AIR IN WALL FOR AIR CABINET UNIT BY OTHERS

ELECTRICAL IN BACKSHELF DETAIL

1. 10120 SINGLE DUPLEX ELECTRIC OUTLET, 20AMP
2. SINGLE JUNCTION BOX WITH BLANK CHROME FINISHED COVER
3. 10220GFWUSB - DOUBLE DUPLEX OUTLETS WITH (2) DUPLEX AND (4) USB PORT 20AMP

SHEET  
NUMBER

02

LOCATION: Naval Branch Health Clinic Building 1600  
Indian Head, MD 20640

DEALER: PATTERSON DENTAL  
Mark Beaupre  
317-413-4556  
mark.beaupre@pattersondental.com


CAD FILE: 218006 Oracle Number: DESIGNER: AJ



5000 SOUTH HALSTED STREET CHICAGO, IL 60609  
PHONE (773) 624-4330 FAX (773) 624-7337



INDIAN HEAD DENTAL LAB, INDIAN HEAD, MD

REVISIONS	
DATE	DISCRIPTION
00/00/2000	
A	
B	
C	
D	
E	

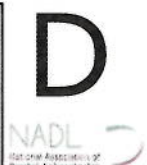








ROOM DIMENSION AND CEILING HEIGHT TO BE VERIFIED BY OTHERS.  
ELECTRICAL, GAS, AIR, AND UTILITIES ABOVE BACKSLASHES IN WALL AS NEEDED, BY CONTRACTOR

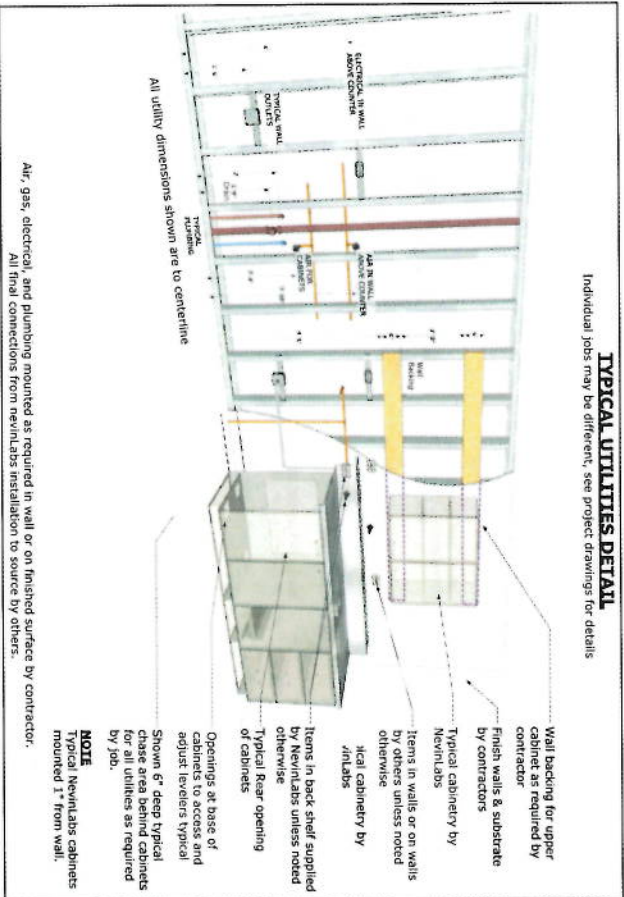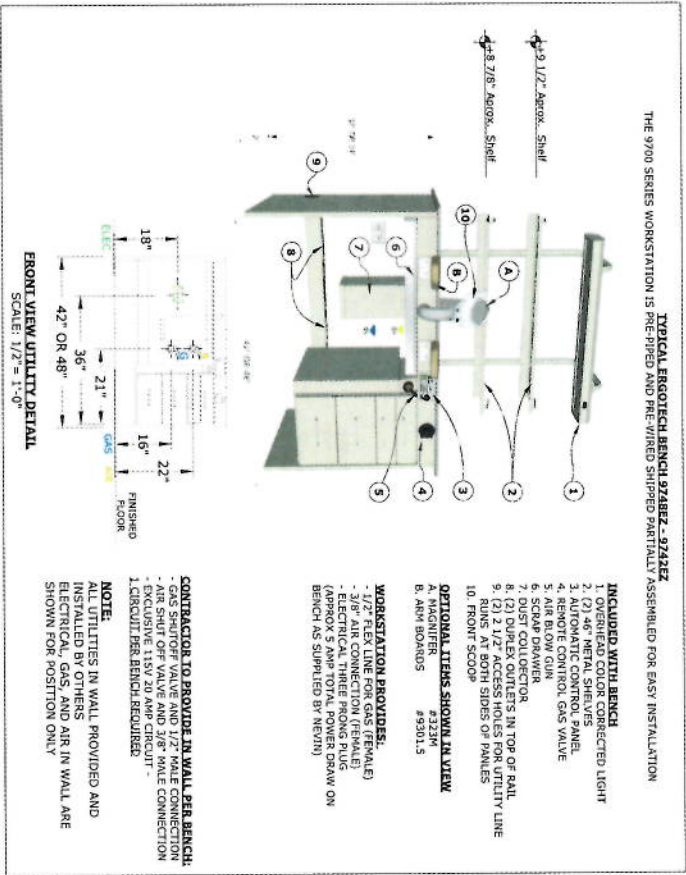

- NOTES:**
- SIZE AND QUANTITIES OF JUNCTION BOXES AND OUTLET ARE SUGGESTED
  - ACTUAL SIZE AND QUANTITIES NEEDED IN FIELD DETERMINED BY OTHERS
  - LOCATIONS OF JUNCTION BOXES AS SHOWN ARE TO CENTERLINE AND TOP OF REAR
  - UTILITIES IN WALL BY CONTRACTOR (AIR, GAS, ELECTRIC) CONTRACTOR SUPPLY SHUT OFF VALVES AND EXCLUSIVE ONE 120V 20 AMP CIRCUIT PER SINGLE STATION.
- REQUIRED UTILITIES IN WALL BY CONTRACTOR FOR WORKSTATIONS (9748Z):**
- 120V 20 AMP EXCLUSIVE CIRCUIT PER BENCH
  - AIR SHUT OFF WITH 3/8" FEMALE CONNECTION
  - GAS SHUT OFF WITH 1/2" FEMALE CONNECTION
- A = AIR  
G = GAS
- ELECTRICAL AND AIR IN WALL ABOVE COUNTER TOP BACK SPLASH AS REQUIRED, BY OTHERS

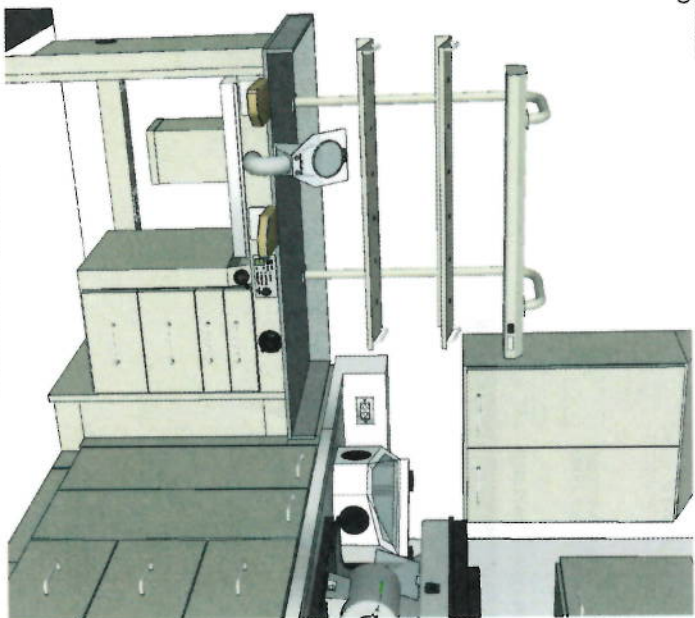

3D VIEW - POSITION 2  
SCALE: N/A

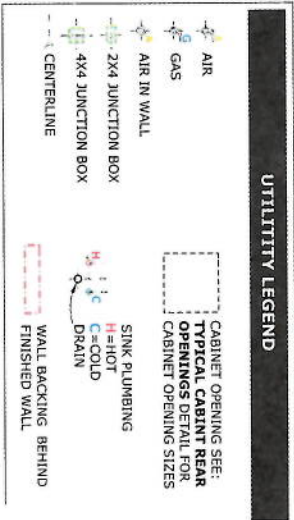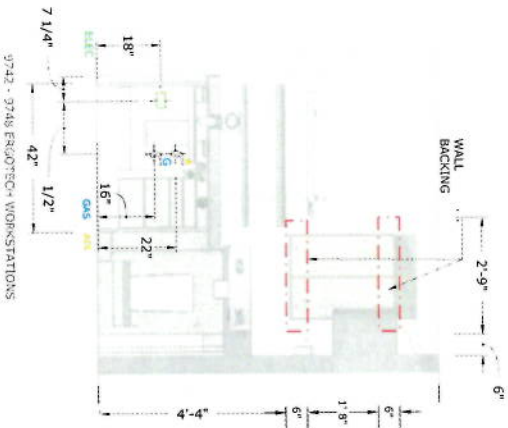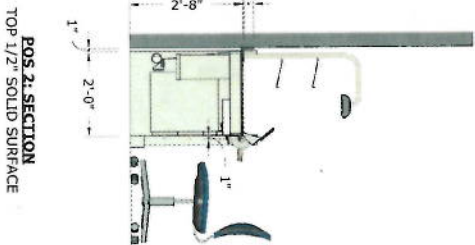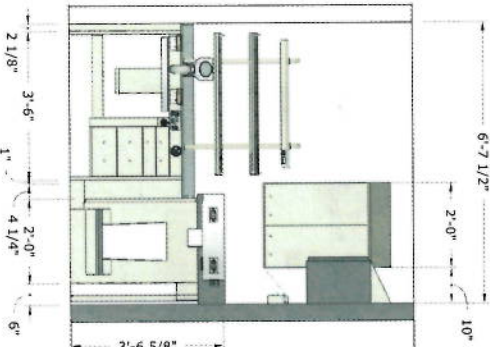

SHEET NUMBER

03

LOCATION: Naval Branch Health Clinic Building 1600  
Indian Head, MD 20640

DEALER: PATTERSON DENTAL  
Mark Beaupre  
317-413-4556  
mark.beaupre@pattersondental.com

CAD FILE: 218006

Oracle Number:

DESIGNER: AJ



5000 SOUTH HALSTED STREET CHICAGO, IL 60609  
PHONE (773) 624-4330 FAX (773) 624-7337

INDIAN HEAD DENTAL LAB, INDIAN HEAD, MD

REVISIONS	
DATE	DISCRIPTION
A 00/00/2000	
B	
C	
D	
E	



NADL  
National Association of  
Dental Laboratories