

checked by TT
on 1/13/15

MEMORANDUM

TO: Mr. Addison Rice
Anderson, Mulholland and Associates

DATE: January 7, 2015

FROM: R. Infante 

FILE: 1412152B

RE: Data Validation
Air samples
SDG: 1412152B

SUMMARY

Full validation was performed on the data for several gas samples analyzed for Methane by method ASTM Method D-1946: Standard Practice for Analysis of Reformed Gas by Gas Chromatography. The samples were collected at the Bristol Myer Squib-Building 5 VI facility, Humacao, PR site on December 09, 2014 and submitted to Eurofins Air Toxics, Inc. of Folsom, California that analyzed and reported the results under delivery group (SDG) 1412152B.

The sample results were assessed according to USEPA data validation guidance documents in the following order of precedence: ASTM Method D-1946: Standard Practice for Analysis of Reformed Gas by Gas Chromatography; Validating Air Samples. Volatile Organic Analysis of Ambient Air in Canisters by Method TO-15, (SOP # HW-31. Revision #4. October, 2006 The QC criteria and data validation actions listed on the data review worksheets are from the primary guidance document, unless otherwise noted.

In general the data is valid as reported and may be used for decision making purposes. The data results are acceptable for use.

SAMPLES

The samples included in the review are listed below

Client Sample ID	Lab. Sample ID	Collected Date	Matrix	Analysis
BSIA-8 (2014)	1412152B-01A	12/09/2014	Air	Methane
BSIA-10 (2014)	1412152B-02A	12/09/2014	Air	Methane
BSIA-6 (2014)	1412152B-03A	12/09/2014	Air	Methane
BSIA-7 (2014)	1412152B-04A	12/09/2014	Air	Methane
BSIA-2 (2014)	1412152B-05A	12/09/2014	Air	Methane
BSIA-1 (2014)	1412152B-04A	12/09/2014	Air	Methane
BSIA-AA (2014)	1412152B-05A	12/09/2014	Air	Methane
BSIA-4 (2014)	1412152B-05A	12/09/2014	Air	Methane

REVIEW ELEMENTS

Sample data were reviewed for the following parameters, where applicable to the method

- Agreement of analysis conducted with chain of custody (COC) form
- Holding time and sample preservation
- Initial and continuing calibrations
- Method blanks/trip blanks/field blank
- Canister cleaning certification criteria
- Field duplicate results
- Laboratory control sample/laboratory control sample duplicate (LCS/LCSD) results
- Quantitation limits and sample results

DISCUSSION

Agreement of Analysis Conducted with COC Request

Sample reports corresponded to the analytical request designated on the chain-of-custody form.

Holding Times and Sample Preservation

Sample preservation was acceptable.

Samples analyzed within method recommended holding time.

Initial and Continuing Calibrations

Methane (ASTM Method D-1946)

The percent relative standard deviations (%RSDs) for response factors (RFs) of all target analytes were within the QC acceptance criteria in the initial calibration. Ongoing accuracy of the instrument was determined by the analysis of a continuing calibration standard.

Method Blank/Trip Blank/Field Blank

Target analytes were not detected in laboratory method blanks for VOCs.

Summa canister met cleaning certification criteria.

Laboratory/Field Duplicate Results

VOCs

Field duplicate were not analyzed as part of this data set. Laboratory duplicate (LCS/LCSD) was used to assess precision. RPD within laboratory/method performance criteria.

LCS/LCSD Results

VOCs

LCS/LCSD (blank spikes) were analyzed by the laboratory associated with this data package. Recoveries and RPD were within laboratory control limits.

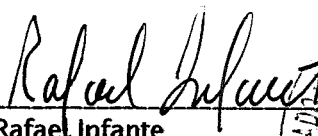
Quantitation Limits and Sample Results


Dilutions were not required with this data set.

Calculations were spot checked.

Certification

The following samples 1412152B-01A; 1412152B-02A; 1412152B-03A; 1412152B-04A; 1412152B-05A; 1412152B-06A; 1412152B-07A; and 1412152B-08A were analyzed following standard procedures accepted by regulatory agencies. The quality control requirements met the methods criteria except in the occasions described in this document. Some of the results were qualified. The results are valid.


Rafael Infante
Chemist License 1888





Air Toxics

Client Sample ID: BSIA- 8 (2014)

Lab ID#: 1412152B-01A

NATURAL GAS ANALYSIS BY MODIFIED ASTM D-1946

File Name:	9121506	Date of Collection: 12/9/14 11:26:00 AM
Dil. Factor:	1.73	Date of Analysis: 12/15/14 01:49 PM

Compound	Rpt. Limit (%)	Amount (%)
Methane	0.00017	0.00033

Container Type: 6 Liter Summa Canister (100% Certified)





Air Toxics

Client Sample ID: BSIA-10 (2014)

Lab ID#: 1412152B-02A

NATURAL GAS ANALYSIS BY MODIFIED ASTM D-1946

File Name:	9121507	Date of Collection:	12/9/14 11:48:00 AM
Dil. Factor:	1.80	Date of Analysis:	12/15/14 02:40 PM

Compound	Rpt. Limit (%)	Amount (%)
Methane	0.00018	0.00032

Container Type: 6 Liter Summa Canister (100% Certified)





Air Toxics

Client Sample ID: BSIA-6 (2014)

Lab ID#: 1412152B-03A

NATURAL GAS ANALYSIS BY MODIFIED ASTM D-1946

File Name:	9121508	Date of Collection: 12/9/14 1:29:00 PM
Dil. Factor:	1.58	Date of Analysis: 12/15/14 03:24 PM

Compound	Rpt. Limit (%)	Amount (%)
Methane	0.00016	0.00090

Container Type: 6 Liter Summa Canister (100% Certified)





Air Toxics

Client Sample ID: BSIA-7 (2014)

Lab ID#: 1412152B-04A

NATURAL GAS ANALYSIS BY MODIFIED ASTM D-1946

File Name:	9121509	Date of Collection:	12/9/14 2:20:00 PM
Dil. Factor:	1.78	Date of Analysis:	12/15/14 03:49 PM

Compound	Rpt. Limit (%)	Amount (%)
Methane	0.00018	0.0012

Container Type: 6 Liter Summa Canister (100% Certified)





Air Toxics

Client Sample ID: BSIA-2 (2014)

Lab ID#: 1412152B-05A

NATURAL GAS ANALYSIS BY MODIFIED ASTM D-1946

File Name:	9121510	Date of Collection:	12/9/14 2:45:00 PM
Dil. Factor:	1.74	Date of Analysis:	12/15/14 04:18 PM

Compound	Rpt. Limit (%)	Amount (%)
Methane	0.00017	0.0010

Container Type: 6 Liter Summa Canister (100% Certified)





Air Toxics

Client Sample ID: BSIA-1 (2014)

Lab ID#: 1412152B-06A

NATURAL GAS ANALYSIS BY MODIFIED ASTM D-1946

File Name:	9121511	Date of Collection:	12/9/14 2:52:00 PM
Dil. Factor:	1.71	Date of Analysis:	12/15/14 04:40 PM
Compound	Rpt. Limit (%)	Amount (%)	
Methane	0.00017	0.0014	

Container Type: 6 Liter Summa Canister (100% Certified)





Air Toxics

Client Sample ID: BSIA-AA (2014)

Lab ID#: 1412152B-07A

NATURAL GAS ANALYSIS BY MODIFIED ASTM D-1946

File Name:	9121512	Date of Collection: 12/9/14 3:08:00 PM
Dil. Factor:	1.67	Date of Analysis: 12/15/14 05:02 PM

Compound	Rpt. Limit (%)	Amount (%)
Methane	0.00017	0.00020

Container Type: 6 Liter Summa Canister (100% Certified)





Air Toxics

Client Sample ID: BSIA-4 (2014)

Lab ID#: 1412152B-08A

NATURAL GAS ANALYSIS BY MODIFIED ASTM D-1946

File Name: 9121513
Dil. Factor: 1.74

Date of Collection: 12/9/14 3:17:00 PM
Date of Analysis: 12/15/14 05:23 PM

Compound	Rpt. Limit (%)	Amount (%)
Methane	0.00017	0.0012

Container Type: 6 Liter Summa Canister (100% Certified)



Sample Transportation Notice

Relinquishing signature on this document indicates that sample is being shipped in compliance with all applicable local, State, Federal, national, and international laws, regulations and ordinances of any kind. Air Toxics Limited assumes no liability with respect to the collection, handling or shipping of these samples. Relinquishing signature also indicates agreement to hold harmless, defend, and indemnify Air Toxics Limited against any claim, demand, or action, of any kind, related to the collection, handling, or shipping of samples. D.O.T. Hotline (800) 467-4922

180 BLUE RAVINE ROAD, SUITE B
FOLSOM, CA 95630-4719
(916) 985-1000 FAX (916) 985-1020

Page 1 of 1

0118 of 0126

Project Manager Terry Taylor
Collected by: (Print and Sign) AK
Company AMAI Email _____
Address 110 Corporate Pk City White Plains State NY Zip 10604
Phone 914-391-1225 Fax _____

Project Info:		Turn Around Time:	Lab Use Only
P.O. # _____	Project # <u>Building S VI</u>	<input checked="" type="checkbox"/> Normal	Pressurized by: _____
Project Name <u>BMS Humacao</u>		<input type="checkbox"/> Rush	Date: _____
specify _____		Pressurization Gas: _____ N ₂ He	

Lab I.D.	Field Sample I.D. (Location)	Can #	Date of Collection	Time of Collection	Analyses Requested	Canister Pressure/Vacuum			
						Initial	Final	Receipt	Final (psi)
01A	BSIA-8(2014)	33654	12/9/14	1126	see notes	28.5	6.5		
02A	BSIA-10(2014)	25261	12/9/14	1148	"	30+	9.5		
03A	BSIA-6(2014)	940	12/9/14	1329	"	27.5	5		
04A	BSIA-7(2014)	33889	12/9/14	1420	"	30	8		
05A	BSIA-2(2014)	3737	12/9/14	1445	"	30+	8		
06A	BSIA-1(2014)	12019	12/9/14	1452	"	30	8		
07A	BSIA-AA(2014)	5695	12/9/14	1508	"	29.5	5		
08A	BSIA-4(2014)	25248	12/9/14	1517	"	30+	7"		

Relinquished by: (signature) <u>[Signature]</u> Date/Time <u>12/9/14 1520</u>	Received by: (signature) <u>[Signature]</u> Date/Time <u>12/10/14 1010</u>	Notes: Acetone, Benzene, Ethylbenzene, toluene, xylene, Methanol, MIBK, Isopropyl Alcohol, Alcohol via TO-15, Methane via ASTM D-1946
Relinquished by: (signature) _____ Date/Time _____	Received by: (signature) _____ Date/Time _____	
Relinquished by: (signature) _____ Date/Time _____	Received by: (signature) _____ Date/Time _____	

Lab Use Only	Shipper Name	Air Bill #	Temp (°C)	Condition	Custody Seals Intact?	Work Order #
	<u>FedEx</u>	<u>77214722 6347</u>	<u>NA</u>	<u>Good</u>	Yes No <u>None</u>	<u>1412152</u>

DATA REVIEW WORKSHEETS

Project Number: 1412152B

Date: 12/09/2014

REVIEW OF VOLATILE ORGANIC PACKAGE

The following guidelines for evaluating volatile organics were created to delineate required validation actions. This document will assist the reviewer in using professional judgment to make more informed decision and in better serving the needs of the data users. The sample results were assessed according to USEPA data validation guidance documents in the following order of precedence: QC criteria from ASTM D-1946 method for measuring permanent gases and light hydrocarbons in refinery and other sources samples using gas chromatography (GC) and a thermal conductivity detector (TCD) and/or flame ionization detection (FID). Validating Air Samples. Volatile Organic Analysis of Ambient Air in Canisters by Method TO-15, (SOP # HW-31. Revision #4. October, 2006). The QC criteria and data validation actions listed on the data review worksheets are from the primary guidance document, unless otherwise noted.

The hardcopied (laboratory name) Eurofins data package received has been reviewed and the quality control and performance data summarized. The data review for VOCs included:

Lab. Project/SDG No.: 1412152B

Sample matrix: Air

No. of Samples: 8

Trip blank No.: -

Field blank No.: -

Equipment blank No.: -

Field duplicate No.: -

☒ Data Completeness

☒ Holding Times

☐ N/A GC/MS Tuning

☐ N/A Internal Standard Performance

☒ Blanks

☐ N/A Surrogate Recoveries

☐ N/A Matrix Spike/Matrix Spike Duplicate

☒ Laboratory Control Spikes

☒ Field Duplicates

☒ Calibrations

☒ Compound Identifications

☒ Compound Quantitation

☒ Quantitation Limits

Overall Comments: Methane by ASTM method D-1946 (modified)

Definition of Qualifiers:

J- Estimated results

U- Compound not detected

R- Rejected data

UJ- Estimated nondetect

Reviewer: Rafael Defaut

Date: 01/07/2015

DATA REVIEW WORKSHEETS

DATA COMPLETENESS

MISSING INFORMATION

DATE LAB. CONTACTED

DATE RECEIVED

This image shows a single sheet of white paper with horizontal blue or grey ruling lines. The lines are evenly spaced and run across the width of the page. There is no handwriting or other markings on the paper.

DATA REVIEW WORKSHEETS

All criteria were met X
Criteria were not met
and/or see below

HOLDING TIMES

The objective of this parameter is to ascertain the validity of the results based on the holding time of the sample from time of collection to the time of analysis.

Complete table for all samples and note the analysis and/or preservation not within criteria

SAMPLE ID	DATE SAMPLED	DATE ANALYZED	pH	ACTION
All samples analyzed within the recommended method holding time				

Criteria

Aqueous samples – 14 days from sample collection for preserved samples ($\text{pH} \leq 2$, 4°C), no air bubbles.

Aqueous samples – 7 days from sample collection for unpreserved samples, 4°C , no air bubbles.

Soil samples- 7 days from sample collection.

Cooler temperature (Criteria: $4 \pm 2^{\circ}\text{C}$): N/A – summa canisters

Actions

If the VOCs vial(s) have air bubbles, estimate positive results (J) and reject nondetects (R).

If the % solids of soil samples is 10-50%, estimate positive results (J) and nondetects (UJ)

If the % solid of soil samples is $< 10\%$, estimate positive results (J) and reject nondetects (R).

If holding times are exceeded but < 14 days beyond criteria, estimate positive results (J) and nondetects (UJ).

If holding times are exceeded but < 28 days beyond criteria, estimate positive results (J) and reject nondetects (R).

If holding times are grossly exceeded (> 28 days beyond criteria), reject all results (R).

If samples were not iced or if the ice were melted ($> 10^{\circ}\text{C}$), estimate positive results (J) and nondetects (UJ).

DATA REVIEW WORKSHEETS

All criteria were met __N/A__
Criteria were not met see below _____

GC/MS TUNING

The assessment of the tuning results is to determine if the sample instrumentation is within the standard tuning QC limits

__N/A__ The BFB performance results were reviewed and found to be within the specified criteria.

__N/A__ BFB tuning was performed for every 24 hours of sample analysis.

If no, use professional judgment to determine whether the associated data should be accepted, qualified or rejected.

List _____ the _____ samples _____ affected:

If mass calibration is in error, all associated data are rejected.

Note: Samples analyzed using GC with either TCD or FID detection.

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below _____

CALIBRATION VERIFICATION

Compliance requirements for satisfactory instrument calibration are established to ensure that the instrument is capable of producing and maintaining acceptable quantitative data.

Date of initial calibration: 05/08/14
 Dates of continuing calibration: 12/15/14
 Instrument ID numbers: GC-9
 Matrix/Level: Air/low

DATE	LAB ID#	FILE	CRITERIA OUT RFs, %RSD, %D, r	COMPOUND	SAMPLES AFFECTED
Initial and continuing calibrations meet method specific requirements. Initial calibration retention times meet method specific requirements.					

Criteria

All RFs must be > 0.05 regardless of method requirements for SPCC.
 All %RSD must be $\leq 15\%$ regardless of method requirements for CCC.
 All %Ds must be $\leq 30\%$ regardless of method requirements for CCC.
 Method TO-15 does not specify criterion for the curve correlation coefficient (r). A limit for r of ≥ 0.995 has therefore been utilized as professional judgment.

Actions

If any compound has an initial RF or a continuing RF of < 0.05 , estimate positive results (J) and reject nondetects (R), regardless of method requirements.
 If any compound has a %RSD $> 15\%$, estimate positive results (J) and use professional judgment to qualify nondetects.
 If any compound has a %RSD $> 90\%$, estimate positive results (J) and reject nondetects (R).
 If any compound has a % D $> 30\%$, estimate positive results (J) and reject nondetects (R).
 If any compound has a % D $> 30\%$, estimate positive results (J) and nondetects (UJ).
 If any compound has a % D $> 90\%$, estimate positive results (J) and reject nondetects (R).
 If any compound has r > 0.995 , estimate positive results and nondetects.

A separate worksheet should be filled for each initial curve

DATA REVIEW WORKSHEETS

All criteria were met X
Criteria were not met
and/or see below

V A. BLANK ANALYSIS RESULTS (Sections 1 & 2)

The assessment of the blank analysis results is to determine the existence and magnitude of contamination problems. The criteria for evaluation of blanks apply only to blanks associated with the samples, including trip, equipment, and laboratory blanks. If problems with any blanks exist, all data associated with the case must be carefully evaluated to determine whether or not there is an inherent variability in the data for the case, or if the problem is an isolated occurrence not affecting other data.

List the contamination in the blanks below. High and low levels blanks must be treated separately.

Laboratory blanks

DATE ANALYZED	LAB ID	LEVEL/MATRIX	COMPOUND	CONCENTRATION UNITS
			All_method_blank_meeth_method_specific_criteria	
			Summa_canisters_met_cleaning_certification_criteria	

Field/Equipment/Trip blank

[illegible]

DATA REVIEW WORKSHEETS

All criteria were met ☒
 Criteria were not met ☐
 and/or see below ☐

V B. BLANK ANALYSIS RESULTS (Section 3)

Blank Actions

Action Levels (ALs) should be based upon the highest concentration of contaminant determined in any blank. Do not qualify any blank with another blank. The ALs for samples which have been diluted should be corrected for the sample dilution factor and/or % moisture, where applicable. No positive sample results should be reported unless the concentration of the compound in the samples exceeds the ALs:

ALs = 10x the amount of common contaminants (methylene chloride, acetone, 2-butanone, and toluene)

ALs = 5x for any other compounds

Specific actions are as follows:

If the concentration is < sample quantitation limit (SQL) and \leq AL, report the compound as not detected (U) at the SQL.

If the concentration is \geq SQL but \leq AL, report the compound as not detected (U) at the reported concentration.

If the concentration is \geq SQL and $>$ AL, report the concentration unqualified.

Notes:

High and low level blanks must be treated separately

Compounds qualified "U" for blank contamination are still considered "hits" when qualifying for calibration criteria.

CONTAMINATION SOURCE/LEVEL	COMPOUND	CONC/UNITS	AL/UNITS	SQL	AFFECTED SAMPLES

DATA REVIEW WORKSHEETS

All criteria were met N/A
 Criteria were not met
 and/or see below

SURROGATE SPIKE RECOVERIES

Laboratory performance of individual samples is established by evaluation of surrogate spike recoveries. All samples are spiked with surrogate compounds prior to sample analysis. The accuracy of the analysis is measured by the surrogate percent recovery. Since the effects of the sample matrix are frequently outside the control of the laboratory and may present relatively unique problems, the validation of data is frequently subjective and demands analytical experience and professional judgment.

List the percent recoveries (%Rs) which do not meet the criteria for surrogate recovery.

Matrix: solid/aqueous

SAMPLE ID	SURROGATE COMPOUND	ACTION
-----------	--------------------	--------

 Surrogate standards not required by the method.

QC Limits* (Air)

_____ LL to UL _____ to _____ to _____ to _____

- * QC limits are laboratory in-house performance criteria, LL = lower limit, UL = upper limit.
- * If QC limits are not available, use limits of 80 – 120 % for aqueous and 70 – 130 % for solid samples.

Actions:

QUALITY	%R < 10%	%R = 10% - LL	%R > UL
Positive results	J	J	J
Nondetects results	R	UJ	Accept

Surrogate action should be applied:

If one or more surrogate in the VOC fraction is out of specification, but has a recovery of > 10%.

If any one surrogate in a fraction shows < 10 % recovery.

DATA REVIEW WORKSHEETS

All criteria were met _____
 Criteria were not met _____
 and/or see below N/A

VII. A MATRIX SPIKE/MATRIX SPIKE DUPLICATE (MS/MSD)

This data is generated to determine long term precision and accuracy in the analytical method for various matrices. This data alone cannot be used to evaluate the precision and accuracy of individual samples. If any % R in the MS or MSD falls outside the designated range, the reviewer should determine if there are matrix effects, i.e. LCS data are within the QC limits but MS/MSD data are outside QC limit.

1. MS/MSD Recoveries and Precision Criteria

The laboratory should use one MS and a duplicate analysis of an unspiked field sample if target analytes are expected in the sample. If target analytes are not expected, MS/MSD should be analyzed.

List the %Rs, RPD of the compounds which do not meet the criteria.

Sample ID: _____ - _____ Matrix/Level: _____ - _____

MS OR MSD	COMPOUND	% R	RPD	QC LIMITS	ACTION
-----------	----------	-----	-----	-----------	--------

MS/MSD are not required as part of ASTM-method D-1946; blank spike used to assess accuracy

- * QC limits are laboratory in-house performance criteria, LL = lower limit, UL = upper limit.
- * If QC limits are not available, use limits of 70 – 130 %.

Actions:

QUALITY	%R < LL	%R > UL
Positive results	J	J
Nondetects results	R	Accept

MS/MSD criteria apply only to the unspiked sample, its dilutions, and the associated MS/MSD samples:

If the % R for the affected compounds were < LL (or 70 %), qualify positive results (J) and nondetects (UJ).

If the % R for the affected compounds were > UL (or 130 %), only qualify positive results (J).

If 25 % or more of all MS/MSD %R were < LL (or 70 %) or if two or more MS/MSD %Rs were < 10%, qualify all positive results (J) and reject nondetects (R).

A separate worksheet should be used for each MS/MSD pair.

DATA REVIEW WORKSHEETS

All criteria were met _____
Criteria were not met _____
and/or see below _____ N/A _____

VII. B MATRIX SPIKE/MATRIX SPIKE DUPLICATE

MS/MSD – Unspiked Compounds

It should be noted that Method TO-15 does not specify a MS/MSD criteria for the unspiked compounds in the sample. A %RSD of < 50% has therefore been utilized as professional judgment.

If all target analytes were spiked in the MS/MSD, this review element is not applicable.

List the %RSD of the compounds which do not meet the criteria.

Sample ID: _____ Matrix/Level/Unit: _____

[illegible]

Actions:

- * If the % RSD > 50, qualify the positive result in the unspiked samples as estimated (J).
* If the % RSD is not calculated (NC) due to nondetected value, use professional judgment to qualify the data.

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

VIII. LABORATORY CONTROL SAMPLE (LCS) ANALYSIS

This data is generated to determine accuracy of the analytical method for various matrices.

1. LCS Recoveries Criteria

Where LCS spiked with the same analyte at the same concentrations as the MS/MSD?
 Yes or No. If no make note in data review memo.

List the %R of compounds which do not meet the criteria

LCS ID	COMPOUND	% R	QC LIMIT
____ LCS/LCSD (Blank spike) analyzed in this data package; recoveries and RPD _____			
____ within laboratory control limits. _____			

- * QC limits are laboratory in-house performance criteria, LL = lower limit, UL = upper limit.
- * If QC limits are not available, use limits of 70 – 130 %.

Actions:

QUALITY	%R < LL	%R > UL
Positive results	J	J
Nondetects results	R	Accept

All analytes in the associated sample results are qualified for the following criteria.

If 25 % of the LCS recoveries were < LL (or 70 %), qualify all positive results (J) and reject nondetects (R).

If two or more LCS were below 10 %, qualify all positive results as (J) and reject nondetects (R).

2. Frequency Criteria:

Where LCS analyzed at the required frequency and for each matrix? Yes or No.

If no, the data may be affected. Use professional judgment to determine the severity of the effect and qualify data accordingly. Discuss any actions below and list the samples affected.

DATA REVIEW WORKSHEETS

All criteria were met _____
 Criteria were not met _____
 and/or see below X

IX. LABORATORY DUPLICATE PRECISION

Sample IDs: _____ - _____

Matrix: _____ - _____

Field/laboratory duplicate samples may be taken and analyzed as an indication of overall precision. These analyses measure both field and lab precision; therefore, the results may have more variability than laboratory duplicates which only laboratory performance. It is also expected that soil duplicate results will have a greater variance than water matrices due to difficulties associated with collecting identical field duplicate samples.

The project QAPP should be reviewed for project-specific information.

Suggested criteria: RPD \pm 25% for air samples. If both samples and duplicate are <5 SQL, the RPD criteria is doubled.

COMPOUND	SQL	SAMPLE CONC.	DUPLICATE CONC.	RPD	ACTION
Field duplicates not analyzed with this data package. Laboratory duplicate (LCS/LCSD) used to assess precision; RPD within laboratory control limits.					

Actions:

Qualify as estimated positive results (J) and nondetects (UJ) for the compound that exceeded the above criteria. For organics, only the sample and duplicate will be qualified.

If an RPD cannot be calculated because one or both of the sample results is not detected, the following actions apply:

If one sample result is not detected and the other is greater than 5x the SQL qualify (J/UJ).

If one sample value is not detected and the other is greater than 5x the SQL and the SQLs for the sample and duplicate are significantly different, use professional judgment to determine if qualification is appropriate.

If one sample value is not detected and the other is less than 5x, use professional judgment to determine if qualification is appropriate.

If both sample and duplicate results are not detected, no action is needed.

DATA REVIEW WORKSHEETS

All criteria were met N/A
 Criteria were not met
 and/or see below

X. INTERNAL STANDARD PERFORMANCE

The assessment of the internal standard (IS) parameter is used to assist the data reviewer in determining the condition of the analytical instrumentation.

List the internal standard area of samples which do not meet the criteria.

- * Area of +40% or -40% of the IS area in the associated calibration standard.
- * Retention time (RT) within ± 0.06 seconds of the IS area in the associated calibration standard.

DATE	SAMPLE ID	IS OUT	IS AREA	ACCEPTABLE RANGE	ACTION
------	-----------	--------	---------	------------------	--------

Internal standard not required by the method. Samples quantified by external standard
method

Actions:

1. IS actions should be applied to the compound quantitated with the out-of-control ISs

QUALITY	IS AREA < -40%		IS AREA > + 40%
Positive results	J		J
Nondetected results	R		ACCEPT

2. If a IS retention time varies more than 0.330 seconds, the chromatographic profile for that sample must be examined to determine if any false positive or negative exists. For shifts of a large magnitude, the reviewer may consider partial or total rejection of the data for the sample fraction.

All criteria were met X

DATA REVIEW WORKSHEETS

Criteria were not met
and/or see below _____

XII. SAMPLE QUANTITATION

The sample quantitation evaluation is to verify laboratory quantitation results. In the space below, please show a minimum of one sample calculation:

1412152B-01A

Methane

RF = 159806712

$$[] = (30278)/(159806712)$$

$$= 0.000189 \% \text{ OK}$$

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

XII. QUANTITATION LIMITS

A. Dilution performed

SAMPLE ID	DILUTION FACTOR	REASONS FOR DILUTION
All samples diluted by a factor of 2		

B. Percent Solids

List samples which have ≤ 50 % solids

Actions:

If the % solids of a soil sample is 10-50%, estimate positive results (J) and nondetects (UJ)

If the % solids of a soil sample is < 10%, estimate positive results (J) and reject nondetects (R)