


File

RECEIVED

JUL 08 2014

AWMD/WRAP-MIRP

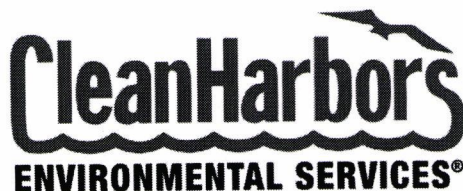

**Memorandum**

**Subject:** Addendum to Phase IV RCRA Facility Investigation Quality Assurance Project Plan for Clean Harbors Kansas (Wichita), Inc., Wichita Kansas

**From:** Martin Smith  
Director of Corrective Measures  
Clean Harbors Environmental Services

**To:** Christine Jump  
EPA Project Manager  
AWMD/RCAP

Clean Harbors is submitting this addendum to the subject report that was approved by EPA on September 9, 2013. The reason for this addendum is Clean Harbors was required to achieve a lower laboratory reporting limit for the herbicide dalapon than was achievable by the laboratory specified in the approved QAPP. Specifically, EPA has requested that select confirmation samples be analyzed for dalapon to demonstrate that soil concentrations in confirmation soil samples are below a Kansas Department of Health and Environment screening level of 0.929 mg/kg. Because this detection limit is not achievable by USEPA Method 8151 (the only dalapon analytical method that Accutest Florida is certified to perform), confirmation samples for dalapon collected as a part of the interim remedial measure at the subject site will be submitted to Test America Seattle where they will be analyzed for dalapon by USEPA Method 8151MS. This method is capable of achieving reporting limits lower than the dalapon interim action objective of 0.929 mg/kg.

Submitted by:

**Martin L Smith**

Digitally signed by Martin L Smith  
DN: cn=Martin L Smith, o=Clean Harbors  
Environmental Services, Inc., ou=Director, Corrective  
Actions and Discontinued Operations,  
email=smith.martin@cleanharbors.com, c=US  
Date: 2014.07.07 19:03:45 -05'00'

Martin Smith  
Director of Corrective Measures  
Clean Harbors Environmental Services

Approved by:

A handwritten signature in blue ink, appearing to read "Christine Jump".

7/8/14

Christine Jump  
EPA Project Manager  
AWMD/RCAP

RCRA



533809

Document Uncontrolled When Printed



Document No. TA-QAM  
Revision No. 7  
Effective Date: 4/23/2014  
Page 1 of 165

## Quality Assurance Manual

TestAmerica Seattle  
5755 8<sup>th</sup> Avenue East  
Tacoma, WA 98424  
Phone No. (253) 922-2310  
Fax No. (253) 922-5047  
www.testamericainc.com

### Copyright Information:

This documentation has been prepared by TestAmerica Laboratories, Inc. and its affiliates ("TestAmerica"), solely for their own use and the use of their customers in evaluating their qualifications and capabilities in connection with a particular project. The user of this document agrees by its acceptance to return it to TestAmerica upon request and not to reproduce, copy, lend, or otherwise disclose its contents, directly or indirectly, and not to use it for any other purpose other than that for which it was specifically provided. The user also agrees that where consultants or other outside parties are involved in the evaluation process, access to these documents shall not be given to said parties unless those parties also specifically agree to these conditions.

THIS DOCUMENT CONTAINS VALUABLE CONFIDENTIAL AND PROPRIETARY INFORMATION. DISCLOSURE, USE OR REPRODUCTION OF THESE MATERIALS WITHOUT THE WRITTEN AUTHORIZATION OF TESTAMERICA IS STRICTLY PROHIBITED. THIS UNPUBLISHED WORK BY TESTAMERICA IS PROTECTED BY STATE AND FEDERAL LAW OF THE UNITED STATES. IF PUBLICATION OF THIS WORK SHOULD OCCUR THE FOLLOWING NOTICE SHALL APPLY:

©COPYRIGHT 2012 TESTAMERICA LABORATORIES, INC. ALL RIGHTS RESERVED.

Facility Distribution No. Uncontrolled

Distributed To: \_\_\_\_\_

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 7  
Effective Date: 4/23/2014  
Page 2 of 165

### Title Page:

## Quality Assurance Manual Approval Signatures

Signatures on File	
Laboratory Director – Heather Prater	Date _____
Quality Assurance Manager – Terri Torres	Date _____
Environmental Health & Safety – Sandra Langway	Date _____
Project Manager Supervisor – Kris Allen	Date _____
Volatile Organic Technical Manager – Bisrat Tadesse	Date _____
Inorganic Technical Manager – Stan Palmquist	Date _____
Semivolatile Technical Manager – Colin McKean	Date _____

Company Confidential & Proprietary



Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 7  
Effective Date: 4/23/2014  
Page 3 of 165

## SECTION 2. TABLE OF CONTENTS

Sec. No.	Title	2009 TNI Standard Reference	ISO/IEC 17025:2005 (E) Reference	Page No.
-	QUALITY ASSURANCE MANUAL	V1M2 Sec. 4.2.8.3		1
1.0	TITLE PAGE			2
2.0	TABLE OF CONTENTS	V1M2 Secs. 4.2.8.3-4.2.8.4		3
3.0	INTRODUCTION, SCOPE AND APPLICABILITY	V1M2 Sec. 4.2.8.4		13
3.1	Introduction And Compliance References	V1M2 Secs. 1.1; 1.2; 2.0; 3.2; 4.1.2; 4.2.4	4.1.2; 4.2.4	13
3.2	Terms And Definitions	V1M2 Secs. 3.0; 4.2.4	4.2.4	13
3.3	Scope / Fields Of Testing	V1M2 Secs. 1.2; 4.2.4	4.1.2; 4.2.4	14
3.4	Management Of The Manual	V1M2 Secs. 4.2.1; 4.2.7; 4.3.2; 4.3.3.3	4.2.1; 4.2.7; 4.3.2; 4.3.3.3	14
4.0	MANAGEMENT REQUIREMENTS	V1M2 Sec. 4		15
4.1	Overview	V1M2 Secs. 4.1.1; 4.1.3; 4.1.5	4.1.1; 4.1.3; 4.1.5; 4.2.6	15
4.2	Roles And Responsibilities	V1M2 Secs. 4.1.4; 4.1.5; 4.1.6; 4.2.1; 4.2.6; 5.2.4	4.1.3; 4.1.5; 4.1.6; 4.2.1; 4.2.6; 5.2.4	15
4.3	Deputies	V1M2 Secs. 4.1.5; 4.1.7.2; 4.2.7	4.1.5; 4.2.7	23
5.0	QUALITY SYSTEM			27
5.1	Quality Policy Statement	V1M2 Secs. 4.1.5; 4.2.2; 4.2.3; 4.2.8.3	4.1.5; 4.2.2; 4.2.3	27
5.2	Ethics And Data Integrity	V1M2 Secs. 4.1.5; 4.1.6; 4.2.2; 4.2.8.1; 5.2.7	4.1.5; 4.2.2	27
5.3	Quality System Documentation	V1M2 Secs. 4.1.5; 4.2.2; 4.2.5	4.2.2; 4.2.5	28
5.4	QA/QC Objectives For The Measurement Of Data	V1M2 Sec. 4.2.2	4.1.5; 4.2.2	29
5.5	Criteria For Quality Indicators			31
5.6	Statistical Quality Control			31
5.7	Quality System Metrics			31
6.0	DOCUMENT CONTROL	V1M2 Secs. 4.2.7; 4.3.1; 4.3.2.2; 4.3.3.3; 4.3.3.4	4.2.7; 4.3.1; 4.3.2.2; 4.3.3.3; 4.3.3.4	32
6.1	Overview			32
6.2	Document Approval And Issue	V1M2 Secs. 4.3.2; 4.3.2.1-4.3.2.3; 4.3.3.1	4.3.2.1; 4.3.2.2; 4.3.2.3; 4.3.3.1	32

Company Confidential & Proprietary

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 7  
Effective Date: 4/23/2014  
Page 4 of 165

Sec. No.	Title	2009 TNI Standard Reference	ISO/IEC 17025:2005 (E) Reference	Page No.
6.3	Procedures For Document Control Policy	V1M2 Secs. 4.3.2.1-4.3.2.2; 4.3.3.1	4.3.2.1; 4.3.2.2; 4.3.3.1	33
6.4	Obsolete Documents	V1M2 Secs. 4.3.2.1-4.3.2.2	4.3.2.1; 4.3.2.2	33
7.0	SERVICE TO THE CLIENT	V1M2 Secs. 4.4.1 - 4.4.4	4.4.1; 4.4.2; 4.4.3; 4.4.4	33
7.1	Overview	V1M2 Secs. 4.4.5; 4.5.5; 5.7.1	4.4.5; 5.7.1	33
7.2	Review Sequence And Key Personnel	V1M2 Sec. 4.4.5	4.4.5	34
7.3	Documentation	V1M2 Sec. 5.7.1	5.7.1	35
7.4	Special Services	V1M2 Secs. 4.7.1-4.7.2	4.7.1; 4.7.2	36
7.5	Client Communication	V1M2 Secs. 4.7.1-4.7.2	4.7.1; 4.7.2	37
7.6	Reporting	V1M2 Secs. 4.7.1-4.7.2	4.7.1; 4.7.2	37
7.7	Client Surveys	V1M2 Secs. 4.7.1-4.7.2	4.7.1; 4.7.2	37
8.0	SUBCONTRACTING OF TESTS	V1M2 Secs. 4.4.3; 4.5.4	4.4.3; 4.5.4	37
8.1	Overview	V1M2 Secs. 4.5.1 - 4.5.3; 4.5.5; 5.3.1	4.5.1; 4.5.2; 4.5.3; 5.3.1	37
8.2	Qualifying And Monitoring Subcontractors	V1M2 Secs. 4.5.1; 4.5.2; 4.5.3; 4.5.5	4.5.1; 4.5.2; 4.5.3	38
8.3	Oversight And Reporting	V1M2 Sec. 4.5.5		39
8.4	Contingency Planning			40
9.0	PURCHASING SERVICES AND SUPPLIES	V1M2 Sec. 4.6.1	4.6.1	42
9.1	Overview	V1M2 Secs. 4.6.2; 4.6.3; 4.6.4	4.6.2; 4.6.3; 4.6.4	42
9.2	Glassware	V1M2 Sec. 5.5.13.1		42
9.3	Reagents, Standards & Supplies	V1M2 Secs. 4.6.2; 4.6.3; 4.6.4	4.6.2; 4.6.3; 4.6.4	42
9.4	Purchase Of Equipment / Instruments / Software			44
9.5	Services			45
9.6	Suppliers			45
10.0	COMPLAINTS	V1M2 Sec. 4.8	4.8	46
10.1	Overview			46
10.2	External Complaints			46
10.3	Internal Complaints			47
10.4	Management Review			47
11.0	CONTROL OF NON-CONFORMING WORK	V1M2 Secs. 4.9.1; 5.10.5	4.9.1; 5.10.5	47
11.1	Overview	V1M2 Secs. 4.9.1; 4.11.3; 4.11.5	4.9.1; 4.11.3; 4.11.5	47

Company Confidential & Proprietary

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 7  
Effective Date: 4/23/2014  
Page 5 of 165

Sec. No.	Title	2009 TNI Standard Reference	ISO/IEC 17025:2005 (E) Reference	Page No.
11.2	Responsibilities And Authorities	V1M2 Secs. 4.9.1; 4.11.3; 4.11.5; 5.2.7	4.9.1; 4.11.3; 4.11.5	48
11.3	Evaluation Of Significance And Actions Taken	V1M2 Secs. 4.9.1; 4.11.3; 4.11.5	4.9.1; 4.11.3; 4.11.5	48
11.4	Prevention Of Nonconforming Work	V1M2 Secs. 4.9.4; 4.11.2	4.9.2; 4.11.2	49
11.5	Method Suspension / Restriction (Stop Work Procedures)	V1M2 Secs. 4.9.1; 4.9.2; 4.11.5	4.9.1; 4.9.2; 4.11.5	49
12.0	CORRECTIVE ACTION	V1M2 Sec. 4.11		50
12.1	Overview	V1M2 Secs. 4.9.2; 4.11.1; 4.11.2	4.9.2; 4.11.1; 4.11.2	50
12.2	General	V1M2 Sec. 4.11.2; 4.11.3	4.11.2; 4.11.3	50
12.3	Closed Loop Corrective Action Process	V1M2 Sec. 4.11.2; 4.11.3; 4.11.4; 4.11.6; 4.11.7; 4.12.2	4.11.2; 4.11.3; 4.11.4; 4.12.2	51
12.4	Technical Corrective Actions	V1M2 Sec. 4.11.6		53
12.5	Basic Corrections	V1M2 Secs. 4.11.1; 4.13.2.3	4.11.1; 4.13.2.3	53
13.0	PREVENTIVE ACTION / IMPROVEMENT	V1M2 Secs. 4.10; 4.12.1; 4.12.2	4.10; 4.12.1; 4.12.2	58
13.1	Overview	V1M2 Secs. 4.15.1; 4.15.2	4.15.1; 4.15.2	58
13.2	Management Of Change			59
14.0	CONTROL OF RECORDS	V1M2 Secs. 4.2.7; 4.13.1.1; 4.13.3	4.2.7; 4.13.1.1	59
14.1	Overview	V1M2 Secs. 4.13.1.1; 4.13.1.2; 4.13.1.3; 4.13.1.4; 4.13.2.1; 4.13.2.2; 4.13.2.3; 4.13.3	4.13.1.1; 4.13.1.2; 4.13.1.3; 4.13.1.4; 4.13.2.1; 4.13.2.2; 4.13.2.3	59
14.2	Technical And Analytical Records	V1M2 Sec. 4.13.2.2 - 4.13.2.3	4.13.2.2; 4.13.2.3	62
14.3	Laboratory Support Activities			63
14.4	Administrative Records			64
14.5	Records Management, Storage And Disposal	V1M2 Sec. 4.13.3		64
15.0	AUDITS			65
15.1	Internal Audits	V1M2 Sec. 4.2.8.1; 4.14; 4.14.1; 4.14.2; 4.14.3; 4.14.5; 5.9.1; 5.9.2	4.14.1; 4.14.2; 4.14.3; 5.9.1; 5.9.2	65
15.2	External Audits	V1M2 Secs. 4.14.2; 4.14.3	4.14.2; 4.14.3; 4.14.4	67

Company Confidential & Proprietary

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 7  
Effective Date: 4/23/2014  
Page 6 of 165

Sec. No.	Title	2009 TNI Standard Reference	ISO/IEC 17025:2005 (E) Reference	Page No.
15.3	Audit Findings	V1M2 Secs. 4.14.2; 4.14.3; 4.14.5		68
16.0	MANAGEMENT REVIEWS	V1M2 Sec. 4.1.6; 4.15; 4.15.1; 4.15.2	4.1.6; 4.15.1; 4.15.2	68
16.1	Quality Assurance Report			68
16.2	Annual Management Review	V1M2 Sec. 4.2.2; 4.15.3	4.2.2	69
16.3	Potential Integrity Related Managerial Reviews			70
17.0	PERSONNEL	V1M2 Secs. 5.2; 5.2.1	5.2.1	70
17.1	Overview	V1M2 Secs. 5.2.2; 5.2.3; 5.2.5	5.2.2; 5.2.3; 5.2.5	70
17.2	Education And Experience Requirements For Technical Personnel	V1M2 Secs. 5.2.1; 5.2.3; 5.2.4	5.2.1; 5.2.3; 5.2.4	71
17.3	Training	V1M2 Sec. 5.2.5	5.2.5	72
17.4	Data Integrity And Ethics Training Program	V1M2 Sec. 4.2.8.1; 5.2.7		73
18.0	ACCOMMODATIONS AND ENVIRONMENTAL CONDITIONS	V1M2 Sec. 5.3		74
18.1	Overview	V1M2 Secs. 5.3.1; 5.3.3; 5.3.4; 5.3.5	5.3.1; 5.3.3; 5.3.4; 5.3.5	74
18.2	Environment	V1M2 Secs. 5.3.1; 5.3.2; 5.3.3; 5.3.4; 5.3.5	5.3.1; 5.3.2; 5.3.3; 5.3.4; 5.3.5	75
18.3	Work Areas	V1M2 Secs. 5.3.3; 5.3.4; 5.3.5	5.3.3; 5.3.4; 5.3.5	75
18.4	Floor Plan			76
18.5	Building Security	V1M2 Sec. 5.3.4	5.3.4	76
19.0	TEST METHODS AND METHOD VALIDATION	V1M2 Sec. 5.4.1	5.4.1	65
19.1	Overview	V1M2 Sec. 5.4.1	5.4.1; 5.4.5.1	76
19.2	Standard Operating Procedures (Sops)	V1M2 Secs. 4.2.8.5; 4.3.3.1; 5.4.2	4.3.3.1; 5.4.2	76
19.3	Laboratory Methods Manual	V1M2 Sec. 4.2.8.5		77

Company Confidential & Proprietary

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 7  
Effective Date: 4/23/2014  
Page 7 of 165

Sec. No.	Title	2009 TNI Standard Reference	ISO/IEC 17025:2005 (E) Reference	Page No.
19.4	Selection Of Methods	V1M2 Secs. 4.13.3; 5.4.1; 5.4.2; 5.4.3; V1M4 Secs. 1.4; 1.5.1; 1.6.1; 1.6.2; 1.6.2.1; 1.6.2.2	5.4.1; 5.4.2; 5.4.3; 5.4.4; 5.4.5.1; 5.4.5.2; 5.4.5.3	77
19.5	Laboratory Developed Methods And Non-Standard Methods	V1M2 Sec. 5.4.2; V1M4 Sec. 1.5.1	5.4.2; 5.4.4; 5.4.5.2; 5.4.5.3	81
19.6	Validation Of Methods	V1M2 Sec. 5.4.2; V1M4 Secs. 1.5.1; 1.5.2; 1.5.2.1; 1.5.2.2; 1.5.3	5.4.2; 5.4.4; 5.4.5.2; 5.4.5.3	81
19.7	Method Detection Limits (mdl) / Limits Of Detection (LOD)	V1M2 Sec. 5.9.3; V1M4 Secs. 1.5.2; 1.5.2.1; 1.5.2.2	5.4.5.3	82
19.8	Instrument Detection Limits (Idl)	V1M2 Sec. 5.9.3		83
19.9	Verification Of Detection And Reporting Limits	V1M2 Sec. 5.9.3; V1M4 Sec. 1.5.2.1		83
19.10	Retention Time Windows	V1M2 Sec. 5.9.3		84
19.11	Evaluation Of Selectivity	V1M2 Sec. 5.9.3; V1M4 Sec. 1.5.4; 1.7.3.6		84
19.12	Estimation Of Uncertainty Of Measurement	V1M2 Sec. 5.1.1; 5.1.2; 5.4.6	5.1.1; 5.1.2; 5.4.6.1; 5.4.6.2; 5.4.6.3	84
19.13	Sample Reanalysis Guidelines	V1M2 Sec 5.9.1	5.9.1	85
19.14	Control Of Data	V1M2 Secs. 5.4.7.1; 5.4.7.2; 5.9.1	5.4.7.1; 5.4.7.2; 5.9.1	85
20.0	EQUIPMENT and CALIBRATIONS	V1M2 Secs. 5.5.4; 5.5.5; 5.5.6	5.5.4; 5.5.5; 5.5.6; 5.6.1	93
20.1	Overview	V1M2 Secs. 5.5.1; 5.5.2; 5.5.3; 5.5.5; 5.5.10	5.5.1; 5.5.2; 5.5.3; 5.5.5; 5.5.10; 5.6.1	93
20.2	Preventive Maintenance	V1M2 Secs. 5.5.1; 5.5.3; 5.5.7; 5.5.9	5.5.1; 5.5.3; 5.5.7; 5.5.9; 5.6.1	93
20.3	Support Equipment	V1M2 Secs. 5.5.10; 5.5.11; 5.5.13.1	5.5.10; 5.5.11; 5.6.2.1.2; 5.6.2.2.1; 5.6.2.2.2	94
20.4	Instrument Calibrations	V1M2 Secs. 5.5.8; 5.5.10; 5.6.3.1; V1M4 Sec. 1.7.1.1; 1.7.2	5.5.8; 5.5.9; 5.5.10; 5.6.1; 5.6.2; 5.6.3.1	96
20.5	Tentatively Identified Compounds (TICS) – GC/MS Analysis			99

Company Confidential & Proprietary

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 7  
Effective Date: 4/23/2014  
Page 8 of 165

Sec. No.	Title	2009 TNI Standard Reference	ISO/IEC 17025:2005 (E) Reference	Page No.
20.6	Gc/Ms Tuning			99
21.0	MEASUREMENT TRACEABILITY			104
21.1	Overview	V1M2 Sec. 5.6.3.1	5.6.2.1.2; 5.6.2.2.2; 5.6.3.1	104
21.2	NIST-Traceable Weights And Thermometers	V1M2 Secs. 5.5.13.1; 5.6.3.1; 5.6.3.2	5.6.3.1; 5.6.3.2	104
21.3	Reference Standards / Materials	V1M2 Secs. 5.6.3.1; 5.6.3.2; 5.6.3.3; 5.6.3.4; 5.6.4.1; 5.6.4.2; 5.9.1; 5.9.3	5.6.3.1; 5.6.3.2; 5.6.3.3; 5.6.3.4; 5.9.1	104
21.4	Documentation And Labeling Of Standards, Reagents, And Reference Materials	V1M2 Secs. 5.6.4.2; 5.9.3		105
22.0	SAMPLING			107
22.1	Overview	V1M2 Secs. 5.7.1; 5.7.3	5.7.1; 5.7.3	107
22.2	Sampling Containers			107
22.3	Definition Of Holding Time			108
22.4	Sampling Containers, Preservation Requirements, Holding Times			108
22.5	Sample Aliquots / Subsampling	V1M2 Sec. 5.7.1	5.7.1	108
23.0	HANDLING OF SAMPLES	V1M2 Sec. 5.8.1	5.8.1	108
23.1	Chain Of Custody (COC)	V1M2 Secs. 5.7.2; 5.7.4; 5.8.4; 5.8.7.5; 5.8.8; 5.9.1	5.7.2; 5.8.4; 5.9.1	108
0	Sample <u>Receipt</u>	V1M2 Secs. 5.8.1; 5.8.2; 5.8.3; 5.8.5; 5.8.7.3; 5.8.7.4; 5.8.7.5	5.8.2; 5.8.3	110
23.3	Sample Acceptance Policy	V1M2 Secs. 5.8.6; 5.8.7.2		111
23.4	Sample Storage	V1M2 Secs. 5.7.4; 5.8.4	5.8.4	111
23.5	Hazardous Samples And Foreign Soils			112
23.6	Sample Shipping	V1M2 Sec. 5.8.2	5.8.2	112
23.7	Sample Disposal			113
24.0	ASSURING THE QUALITY OF TEST RESULTS			119
24.1	Overview	V1M2 Secs. 5.9.2; 5.9.3	5.9.2	119
24.2	Controls	V1M2 Secs. 5.9.2; 5.9.3	5.9.2	119

Company Confidential & Proprietary

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 7  
Effective Date: 4/23/2014  
Page 9 of 165

Sec. No.	Title	2009 TNI Standard Reference	ISO/IEC 17025:2005 (E) Reference	Page No.
24.3	Negative Controls	V1M2 Secs. 5.9.2; 5.9.3 V1M4 Secs. 1.7.3; 1.7.3.1; 1.7.4.1	5.9.2	119
24.4	Positive Controls	V1M2 Secs. 5.9.2; 5.9.3 V1M4 Secs. 1.7.3; 1.7.3.2; 1.7.3.2.1; 1.7.3.2.2; 1.7.3.2.3	5.9.2	120
24.5	Sample Matrix Controls	V1M2 Secs. 5.9.2; 5.9.3 V1M4 Secs. 1.7.3; 1.7.3.3; 1.7.3.3.1; 1.7.3.3.2; 1.7.3.3.3	5.9.2	122
24.6	Acceptance Criteria (Control Limits)	V1M2 Sec. 5.9.3. V1M4 Secs. 1.7.4.2; 1.7.4.3		122
24.7	Additional Procedures To Assure Quality Control	V1M2 Sec. 5.9.3. V1M4 Sec. 1.7.3.4		125
25.0	REPORTING RESULTS			125
25.1	Overview	V1M2 Secs. 5.10.1; 5.10.2; 5.10.8	5.10.1; 5.10.2; 5.10.8	125
25.2	Test Reports	V1M2 Secs. 5.10.1; 5.10.2; 5.10.3.1; 5.10.3.2; 5.10.5; 5.10.6; 5.10.7; 5.10.8; 5.10.10; 5.10.11	5.10.1; 5.10.2; 5.10.3.1; 5.10.3.2; 5.10.5; 5.10.6; 5.10.7; 5.10.8	125
25.3	Reporting Level Or Report Type	V1M2 Secs. 5.10.1; 5.10.7; 5.10.8	5.10.1; 5.10.7; 5.10.8	128
25.4	Supplemental Information For Test	V1M2 Secs. 5.10.1; 5.10.3.1; 5.10.5	5.10.1; 5.10.3.1; 5.10.5	128
25.5	Environmental Testing Obtained From Subcontractors	V1M2 Secs. 4.5.5; 5.10.1; 5.10.6	5.10.1; 5.10.6	129
25.6	Client Confidentiality	V1M2 Secs. 4.1.5; 5.10.7	4.1.5; 5.10.7	129
25.7	Format Of Reports	V1M2 Sec. 5.10.8	5.10.8	130
25.8	Amendments To Test Reports	V1M2 Sec. 5.10.9	5.10.1; 5.10.9	130
25.9	Policies On Client Requests For Amendments	V1M2 Secs. 5.9.1; 5.10.9	5.9.1; 5.10.1; 5.10.5; 5.10.9	130

Company Confidential & Proprietary

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 7  
Effective Date: 4/23/2014  
Page 10 of 165

## LIST OF TABLES

Table No.	Title	2009 TNI Standard Reference	ISO/IEC 17025:2005 (E) Reference	Page No.
12-1	Example – General Corrective Action Procedures	V1M2 Sec. 4.11.6. V1M4 Sec. 1.7.4.1	4.11.2	55
14-1	Record Index		4.13.1.1	59
14-2	Example: Special Record Retention Requirements			61
15-1	Types Of Internal Audits And Frequency		4.14.1	66
20-1	Example: Instrumentation List		5.5.4; 5.5.5	99
20-2	Example: Software List			100
20-3	Example: Schedule Of Routine Maintenance			101
24-1	Example – Negative Controls			119
24-3	Sample Matrix Control			122

## LIST OF FIGURES

Figure No.	Title	2009 TNI Standard Reference	ISO/IEC 17025:2005(E) Reference	Page No.
4-1	Corporate And Laboratory Organization Chart	V1M2 Sec. 4.1.5	4.1.3; 4.1.5; 4.2.6	255
8-1	Example - Subcontracted Sample Form			41
12-1	Example - Non Conformance Memo			54
19-1	Example - Demonstration Of Capability Documentation			91
19-2	Example: Work Flow			92
23-1	Example: Chain Of Custody (COC)			114
23-2	Example: Sample Acceptance Policy	V1M2 Sec. 5.8.6; 5.8.7.1. V1M4 Sec. 1.7.5		115
23-3	Example: Cooler Receipt Form		5.8.3	118

Company Confidential & Proprietary

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 7  
Effective Date: 4/23/2014  
Page 11 of 165

#### LIST OF APPENDICES

Appendix No.	Title	Page No.
1	Laboratory Floor Plan	132
2	Glossary/Acronyms	133
3	Laboratory Certifications, Accreditations, Validations	141
4	Summary Of Calibration, QC And Corrective Action Procedures	142

#### REFERENCED CORPORATE SOPs AND POLICIES

SOP / Policy Reference	Title
CA-Q-S-001	Solvent and Acid Lot Testing and Approval
CA-Q-S-002	Acceptable Manual Integration Practices
CA-Q-S-004	Method Compliance & Data Authenticity Audits
CA-Q-S-006	Detection Limits
CA-Q-S-008	Management Systems Review
CA-Q-S-009	Root Cause Analysis
CW-Q-S-001	Corporate Document Control and Archiving
CW-Q-S-002	Writing a Standard Operating Procedure (SOPs)
CW-L-S-002	Internal Investigation of Potential Data Discrepancies and Determination for Data Recall
CA-L-S-002	Subcontracting Procedures
CW-L-P-004	Ethics Policy
CA-L-P-002	Contract Compliance Policy
CW-F-P-002	Authorization Matrix
CW-F-P-004	Procurement and Contracts Policy
CA-C-S-001	Work Sharing Process
CA-T-P-001	Qualified Products List
CW-F-S-007	Controlled Purchases Policy
CW-F-S-018	Vendor Selection
CA-Q-M-002	Corporate Quality Management Plan
CW-E-M-001	Corporate Environmental Health & Safety Manual

Company Confidential & Proprietary

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 7  
Effective Date: 4/23/2014  
Page 12 of 165

#### REFERENCED LABORATORY SOPs

SOP Reference	Title
TA-QA-0528	Document Control
TA-QA-0529	Complaint Resolution
TA-QA-0530	Management of Change
TA-QA-0506	Archiving Reports and Report File Maintenance
TA-QA-0608	Employee Training Procedures
TA-QA-0500	Standard Operating Procedures
TA-QA-0617	Analyst Demonstration of Capability
TA-QA-0602	Method Detection Limit Studies and Limits of Detection
TA-IP-0226	Multi Incremental Subsampling of Soils and Sediments
TA-QA-0028	Subsampling of Solid Samples
TA-QA-0001	Sample Receipt & Login
TA-QA-0610	Laboratory Corrective Action Procedures
TA-QA-0635	Data Review
TA-QA-0014	Selection and Use of Laboratory Balances
TA-QA-0024	Use, Calibration and Maintenance of Laboratory Thermometers
TA-QA-0016	Volumetric Verification
TA-QA-0610	Laboratory Corrective Action Procedures
TA-QA-0619	Receipt, Storage and Verification of Standards
TA-QA-0600	Quality Control Charting and Establishing

Company Confidential & Proprietary



### SECTION 3. INTRODUCTION, SCOPE AND APPLICABILITY

#### 3.1 Introduction and Compliance References

TestAmerica Seattle's Quality Assurance Manual (QAM) is a document prepared to define the overall policies, organization objectives and functional responsibilities for achieving TestAmerica's data quality goals. The laboratory maintains a local perspective in its scope of services and client relations and maintains a national perspective in terms of quality.

The QAM has been prepared to assure compliance with The NELAC Institute (TNI) Standard, dated 2009, Volume 1 Modules 2 and 4, and ISO/IEC Guide 17025:2005(E). In addition, the policies and procedures outlined in this manual are compliant with TestAmerica's Corporate Quality Management Plan (CQMP) and the various accreditation and certification programs listed in Appendix 3. The CQMP provides a summary of TestAmerica's quality and data integrity system. It contains requirements and general guidelines under which all TestAmerica facilities shall conduct their operations.

The QAM has been prepared to be consistent with the requirements of the following documents:

- EPA 600/4-88/039, Methods for the Determination of Organic Compounds in Drinking Water, EPA, Revised July 1991.
- EPA 600/R-95/131, Methods for the Determination of Organic Compounds in Drinking Water, Supplement III, EPA, August 1995.
- EPA 600/4-79-019, Handbook for Analytical Quality Control in Water and Wastewater Laboratories, EPA, March 1979.
- Test Methods for Evaluating Solid Waste Physical/Chemical Methods (SW846), Third Edition, September 1986, Final Update I, July 1992, Final Update IIA, August 1993, Final Update II, September 1994; Final Update IIB, January 1995; Final Update III, December 1996; Final Update IV, January 2008.
- U.S. Department of Defense, Quality Systems Manual for Environmental Laboratories, Version 4.2, October 2010.
- Federal Register, 40 CFR Parts 136, 141, 172, 173, 178, 179 and 261.
- Manual for the Certification of Laboratories Analyzing Drinking Water (EPA 815-R-05-004, January 2005) (DW labs only)
- APHA, Standard Methods for the Examination of Water and Wastewater, 18<sup>th</sup> Edition, 19<sup>th</sup>, 20<sup>th</sup>, 21<sup>st</sup>, and on-line Editions.
- Marine Protection, Research, and Sanctuaries Act (MPRSA).
- Toxic Substances Control Act (TSCA).

#### 3.2 Terms and Definitions

A Quality Assurance Program is a company-wide system designed to ensure that data produced by the laboratory conforms to the standards set by state and/or federal regulations. The program functions at the management level through company goals and management policies, and at the analytical level through Standard Operating Procedures (SOPs) and quality

control. The TestAmerica program is designed to minimize systematic error, encourage constructive, documented problem solving, and provide a framework for continuous improvement within the organization.

Refer to Appendix 2 for the Glossary/Acronyms.

#### 3.3 Scope / Fields of Testing

The laboratory analyzes a broad range of environmental and industrial samples every month. Sample matrices vary among drinking water, effluent water, groundwater, hazardous waste, sludge, sediments and soils. The Quality Assurance Program contains specific procedures and methods to test samples of differing matrices for chemical and physical parameters. The Program also contains guidelines on maintaining documentation of analytical processes, reviewing results, servicing clients and tracking samples through the laboratory. The technical and service requirements of all analytical requests are thoroughly evaluated before commitments are made to accept the work. Measurements are made using published reference methods or methods developed and validated by the laboratory.

The methods covered by this manual include the most frequently requested methodologies needed to provide analytical services in the United States and its territories. The specific list of test methods used by the laboratory can be found in the Laboratory Information Management System (LIMS). The approach of this manual is to define the minimum level of quality assurance and quality control necessary to meet these requirements. All methods performed by the laboratory shall meet these criteria as appropriate. In some instances, quality assurance project plans (QAPPs), project specific data quality objectives (DQOs) or local regulations may require criteria other than those contained in this manual. In these cases, the laboratory will abide by the requested criteria following review and acceptance of the requirements by the Laboratory Director and the Quality Assurance (QA) Manager. In some cases, QAPPs and DQOs may specify less stringent requirements. The Laboratory Director and the QA Manager must determine if it is in the lab's best interest to follow the less stringent requirements.

#### 3.4 Management of the Manual

##### 3.4.1 Review Process

The template on which this manual is based is reviewed annually by Corporate Quality Management Personnel to assure that it remains in compliance with Section 3.1. This manual itself is reviewed **annually** by senior laboratory management to assure that it reflects current practices and meets the requirements of the laboratory's clients and regulators as well as the CQMP. Occasionally, the manual may need changes in order to meet new or changing regulations and operations. The QA Manager will review the changes in the normal course of business and incorporate changes into revised sections of the document. All updates will be reviewed by the senior laboratory management staff. The laboratory updates and approves such changes according to our Document Control procedures (refer to SOP No. TA-QA-0528).

## SECTION 4. MANAGEMENT REQUIREMENTS

### 4.1 Overview

TestAmerica Seattle is a local operating unit of TestAmerica Laboratories, Inc.. The organizational structure, responsibilities and authorities of the corporate staff of TestAmerica Laboratories, Inc. are presented in the CQMP. The laboratory has day-to-day independent operational authority overseen by corporate officers (e.g., President, Chief Executive Officer, Corporate Quality, etc.). The laboratory operational and support staff work under the direction of the Laboratory Director. The organizational structure for both Corporate & TestAmerica Seattle is presented in Figure 4-1.

### 4.2 Roles and Responsibilities

In order for the Quality Assurance Program to function properly, all members of the staff must clearly understand and meet their individual responsibilities as they relate to the quality program. The following descriptions briefly define each role in its relationship to the Quality Assurance Program.

#### 4.2.1 Additional Requirements for Laboratories

The responsibility for quality resides with every employee of the laboratory. All employees have access to the QAM, are trained to this manual, and are responsible for upholding the standards therein. Each person carries out his/her daily tasks in a manner consistent with the goals and in accordance with the procedures in this manual and the laboratory's SOPs. Role descriptions for Corporate personnel are defined in the CQMP. This manual is specific to the operations of TestAmerica's Seattle laboratory.

#### 4.2.2 Laboratory Director

TestAmerica Seattle's Laboratory Director is responsible for the overall quality, safety, financial, technical, human resource and service performance of the whole laboratory and reports to their respective GM. The Laboratory Director provides the resources necessary to implement and maintain an effective and comprehensive Quality Assurance and Data Integrity Program.

Specific responsibilities include, but are not limited to:

- Provides one or more technical managers for the appropriate fields of testing. The name(s) of the Technical Director will be included in the national database. If the Technical Director is absent for a period of time exceeding 15 consecutive calendar days, the Laboratory Director must designate another full time staff member meeting the qualifications of the Technical Director to temporarily perform this function. If the absence exceeds 65 consecutive calendar days, the primary accrediting authority must be notified in writing.
- Ensures that all analysts and supervisors have the appropriate education and training to properly carry out the duties assigned to them and ensures that this training has been documented.
- Ensures that personnel are free from any commercial, financial and other undue pressures that might adversely affect the quality of their work.

- Ensures TestAmerica's human resource policies are adhered to and maintained.
- Ensures that sufficient numbers of qualified personnel are employed to supervise and perform the work of the laboratory.
- Ensures that appropriate corrective actions are taken to address analyses identified as requiring such actions by internal and external performance or procedural audits. Procedures that do not meet the standards set forth in the QAM or laboratory SOPs may be temporarily suspended by the Laboratory Director.
- Reviews and approves the laboratory specific QAM, policies, SOPs prior to their implementation and ensures all approved procedures are implemented and adhered to.
- Pursues and maintains appropriate laboratory certification and contract approvals. Supports ISO 17025 requirements. Organizes bid activities for prospective new projects and clients.
- Ensures client specific reporting and quality control requirements are met.
- Captains the management team, consisting of the QA Manager, the Technical Director(s), and the Operations Manager as direct reports.
- Annually assesses the effectiveness of the QMP and QAM within the laboratory.

#### 4.2.3 Quality Assurance (QA) Manager or Designee

The QA Manager has responsibility and authority to ensure the continuous implementation of the quality system.

The QA Manager reports directly to the Laboratory Director and has access to Corporate QA for advice and resources. This position is able to evaluate data objectively and perform assessments without outside (e.g., managerial) influence. Corporate QA may be used as a resource in dealing with regulatory requirements, certifications and other quality assurance related items. The QA Manager directs the activities of the QA officers to accomplish specific responsibilities, which include, but are not limited to:

- Compliance with ISO 17025
- Compliance with DoD ELAP
- Serves as the focal point for QA/QC in the laboratory.
- Having functions independent from laboratory operations for which he/she has quality assurance oversight.
- Maintaining and updating the QAM.
- Monitoring and evaluating laboratory certifications; scheduling proficiency testing samples.
- Monitoring and communicating regulatory changes that may affect the laboratory to management.
- Training and advising the laboratory staff on quality assurance/quality control procedures that are pertinent to their daily activities.
- Have documented training and/or experience in QA/QC procedures and the laboratory's Quality System.

## Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 7  
Effective Date: 4/23/2014  
Page 17 of 165

- Having a general knowledge of the analytical test methods for which data audit/review is performed (and/or having the means of getting this information when needed).
- Arranging for or conducting internal audits on quality systems and the technical operation.
- The laboratory QA Manager will maintain records of all ethics-related training, including the type and proof of attendance.
- Maintain, improve, and evaluate the corrective action database and the corrective and preventive action systems.
- Notifying laboratory management of deficiencies in the quality system and ensuring corrective action is taken. Procedures that do not meet the standards set forth in the QAM or laboratory SOPs shall be investigated following procedures outlined in Section 12 and if deemed necessary may be temporarily suspended during the investigation.
- Objectively monitors standards of performance in quality control and quality assurance without outside (e.g., managerial) influence.
- Coordinating of document control of SOPs, MDLs, control limits, and miscellaneous forms and information.
- Review a percentage of all final data reports for internal consistency. Review of Chain of Custody (COC), correspondence with the analytical request, batch QC status, completeness of any corrective action statements, 5% of calculations, format, holding time, sensibility and completeness of the project file contents.
- Review of external audit reports and data validation requests.
- Follow-up with audits to ensure client QAPP requirements are met.
- Establishment of reporting schedule and preparation of various quality reports for the Laboratory Director, clients and/or Corporate QA.
- Development of suggestions and recommendations to improve quality systems.
- Research of current state and federal requirements and guidelines.
- Captains the QA team to enable communication and to distribute duties and responsibilities.
- Ensuring Communication & monitoring standards of performance to ensure that systems are in place to produce the level of quality as defined in this document.
- Evaluation of the thoroughness and effectiveness of training.

### 4.2.4 Technical Managers or Designees

The Technical Managers report directly to the Laboratory Director. He/she is accountable for all analyses and analysts under their experienced supervision. The scope of responsibility ranges from the new-hire process and existing technology through the ongoing training and development programs for existing analysts and new instrumentation. Specific responsibilities include, but are not limited to:

- Compliance with ISO 17025
- Compliance with DoD ELAP

Company Confidential & Proprietary

## Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 7  
Effective Date: 4/23/2014  
Page 18 of 165

- Exercises day-to-day supervision of laboratory operations for the appropriate field of accreditation and reporting of results. Coordinating, writing, and reviewing preparation of all test methods, i.e., SOPs, with regard to quality, integrity, regulatory and optimum and efficient production techniques, and subsequent analyst training and interpretation of the SOPs for implementation and unusual project samples. He/she insures that the SOPs are properly managed and adhered to at the bench. He/she develops standard costing of SOPs to include supplies, labor, overhead, and capacity (design vs. demonstrated versus first-run yield) utilization.
- Reviewing and approving, with input from the QA Manager, proposals from marketing, in accordance with an established procedure for the review of requests and contracts. This procedure addresses the adequate definition of methods to be used for analysis and any limitations, the laboratory's capability and resources, the client's expectations. Differences are resolved before the contract is signed and work begins. A system documenting any significant changes is maintained, as well as pertinent discussions with the client regarding their requirements or the results of the analyses during the performance of the contract. All work subcontracted by the laboratory must be approved by the client. Any deviations from the contract must be disclosed to the client. Once the work has begun, any amendments to the contract must be discussed with the client and so documented.
- Monitoring the validity of the analyses performed and data generated in the laboratory. This activity begins with reviewing and supporting all new business contracts, insuring data quality, analyzing internal and external non-conformances to identify root cause issues and implementing the resulting corrective and preventive actions, facilitating the data review process (training, development, and accountability at the bench), and providing technical and troubleshooting expertise on routine and unusual or complex problems.
- Providing training and development programs to applicable laboratory staff as new hires and, subsequently, on a scheduled basis. Training includes instruction on calculations, instrumentation management to include troubleshooting and preventive maintenance.
- Enhancing efficiency and improving quality through technical advances and improved LIMS utilization. Capital forecasting and instrument life cycle planning for second generation methods and instruments as well as asset inventory management.
- Coordinating sample management from "cradle to grave," insuring that no time is lost in locating samples.
- Scheduling all QA/QC-related requirements for compliance, e.g., MDLs, etc..
- Captains department personnel to communicate quality, technical, personnel, and instrumental issues for a consistent team approach.
- Coordinates audit responses with the QA Manager.

### 4.2.5 LIMS Administration/IT Support

LIMS Administration and IT support are managed at the corporate level. In general, the duties of the corporate LIMS administration team consist of:

- Establishes and maintains the laboratory information system (LIMS) for tracking all samples in the laboratory.
- Updates and enhances LIMS.

Company Confidential & Proprietary

- Develops expertise in the requirements described in Good Automated Laboratory Practices (GALP) EPA 2185, 1995 Edition, in order to ensure compliance.
- Programs and tests software modifications/changes.
- Coordinates testing to ensure that all LIMS software accurately performs its intended functions. Testing is performed and documented after installation or when modifications/changes are made.
- Maintains historical files of software, software operating procedures (manuals), software changes/modifications (Change Log) and software version numbers.
- Maintains log of repairs and service performed on LIMS hardware.
- Develops and verifies security practices to assure the integrity of LIMS data. Identifies threats, potential threats, and future threats.
- Maintains awareness of any environmental conditions of the facility housing the LIMS that may compromise LIMS raw data and informs management.
- LIMS database back-up once daily.

#### 4.2.6 Hazardous Waste Coordinator

The Hazardous Waste Coordinator reports directly to the Laboratory Director. The duties consist of:

- Managing laboratory generated hazardous waste in accordance with appropriate regulations.
- Staying current with the hazardous waste regulations.
- Continuing training on hazardous waste issues.
- Reviewing and updating annually the Hazardous Waste Contingency Plan in the Environmental Health & Safety Manual.
- Auditing the staff with regard to compliance with the Hazardous Waste Contingency Plan.
- Contacting the hazardous waste subcontractors for review of procedures and opportunities

#### 4.2.8 Department Managers

The Department Manager duties at TestAmerica Seattle are performed by the Technical Managers. The responsibilities of the department managers are listed below:

- Ensure that analysts in their department adhere to applicable SOPs and the QA Manual. They perform frequent SOP and QA Manual review to determine if analysts are in compliance and if new, modified, and optimized measures are feasible and should be added to these documents.
- With regard to analysts, participates in the selection, training (as documented in Section 8.1), development of performance objectives and standards of performance, appraisal (measurement of objectives), scheduling, counseling, discipline, and motivation of analysts and documents these activities in accordance with systems developed by the QA and Personnel Departments. They evaluate staffing sufficiency and overtime needs. Training consists of familiarization with SOP, QC, Safety, and computer systems.

- Encourage the development of analysts to become cross-trained in various methods and/or operate multiple instruments efficiently while performing maintenance and documentation, self-supervise, and function as a department team.
- Provide guidance to analysts in resolving problems encountered daily during sample prep/analysis in conjunction with the QA Manager. Each is responsible for 100% of the data review and documentation, non-conformance and CPAR issues, the timely and accurate completion of performance evaluation samples and MDLs, for his/her department.
- Works with the Project Managers and Analysts/Technicians to ensure that the requirements of project are met in a timely manner.
- Ensure all logbooks are maintained, current, and properly labeled or archived.
- Responsible for meeting quality requirements and reporting all non-conformance conditions to the QA Manager and/or Laboratory Director.
- Ensure that preventive maintenance is performed and documented on instrumentation as detailed in the QA Manual or SOPs. He is responsible for developing and implementing a system for preventive maintenance, troubleshooting, and repairing or arranging for repair of instruments.
- Maintain adequate and valid inventory of reagents, standards, spare parts, and other relevant resources required to perform daily analysis.
- Achieve optimum turnaround time on analyses and compliance with holding times.
- Conduct efficiency and cost control evaluations on an ongoing basis to determine optimization of labor, supplies, overtime, first-run yield, capacity (designed vs. demonstrated), second- and third-generation production techniques/instruments, and long-term needs for budgetary planning.
- Develop, implement, and enhance calibration programs.
- Provide written responses and facilitates implementing corrective actions for external and internal audit issues.

#### 4.2.9 Laboratory Analysts

Laboratory analysts are responsible for conducting analysis and performing all tasks assigned to them by the department manager, group leader or supervisor. The responsibilities of the analysts are listed below:

- Perform analyses by adhering to analytical and quality control protocols prescribed by current SOPs, this QA Manual, and project-specific plans honestly, accurately, timely, safely, and in the most cost-effective manner.
- Document standard and sample preparation, instrument calibration and maintenance, data calculations, sample matrix effects, and any observed non-conformance on worklists, benchsheets, lab notebooks and/or the Non-Conformance Database.
- Report all non-conformance situations, instrument problems, matrix problems and QC failures, which might affect the reliability of the data, to their supervisor, the Technical Director, and/or the QA Manager or member of QA staff.
- Ensures sample and data integrity by adhering to internal chain-of-custody procedures.

- Perform 100% review of the data generated prior to entering and submitting for secondary level review.
- Suggest method improvements to their supervisor, the Technical Director, and the QA Manager. These improvements, if approved, will be incorporated. Ideas for the optimum performance of their assigned area, for example, through the proper cleaning and maintenance of the assigned instruments and equipment, are encouraged.
- Work cohesively as a team in their department to achieve the goals of accurate results, optimum turnaround time, cost effectiveness, cleanliness, complete documentation, and personal knowledge of environmental analysis.

#### 4.2.10 Safety Officer

The Safety Officer reports to the Laboratory Director and ensures that systems are maintained for the safe operation of the laboratory. The Safety Officer is responsible to:

- Conduct ongoing, necessary safety training and conduct new employee safety orientation.
- Assist in developing and maintaining the Chemical Hygiene/Safety Manual.
- Administer dispersal of all Material Safety Data Sheet (MSDS) information.
- Perform regular chemical hygiene and housekeeping instruction.
- Give instruction on proper labeling and practice.
- Serve as chairman of the laboratory safety committee.
- Provide and train personnel on protective equipment.
- Oversee the inspection and maintenance of general safety equipment – fire extinguishers, safety showers, eyewash fountains, etc. and ensure prompt repairs as needed.
- Supervise and schedule fire drills and emergency evacuation drills.
- Determine what initial and subsequent exposure monitoring, if necessary to determine potential employee exposure to chemicals used in the laboratory.
- When determined necessary, conduct exposure monitoring assessments.
- Determine when a complaint of possible over-exposure is "reasonable" and should be referred for medical consultation.
- Assist in the internal and external coordination of the medical consultation/monitoring program conducted by TestAmerica's medical consultants.
- Manages facility maintenance.

#### 4.2.11 Sample Receiving Staff

The Sample Receiving Technicians report to the Sample Receiving Supervisor; the Supervisor reports to the Operations Manager. The responsibilities are outlined below:

- Ensures implementation of proper sample receipt procedures, including maintenance of chain-of-custody.
- Reports nonconformances associated with condition-upon-receipt of samples.
- Logs incoming samples into the LIMS.
- Ensures that all samples are stored in the proper environment.
- Ensure the verification of data entry from login.
- Responsible for meeting quality requirements including documenting preservation.
- Responsible for ensuring the timely and correct shipment of sample containers, including proper preservatives and instructions, to clients
- Assists Environmental Health and Safety staff with sample disposal.

#### 4.2.12 Project Manager Supervisor

The Project Manager Supervisor reports to the Client Services Director and serves as the interface between the laboratory's technical departments and the laboratory's clients. The staff consists of the Project Management team. With the overall goal of total client satisfaction, the functions of this position are outlined below:

- Technical training and growth of the Project Management team.
- Technical liaison for the Project Management team.
- Human resource management of the Project Management team.
- Has signature authority for laboratory reports.
- Assesses and assures customer satisfaction.
- Provides feedback to management on changing customer needs.
- Works with the Department Managers and/or Analysts/Technicians to ensure the requirements of projects are met in a timely manner.
- Organizes bid activities for prospective new projects and clients.



#### 4.2.13 Project Managers

The Project Managers report to the Customer Service Manager and serve as the interface between the laboratory and the clients. With the overall goal of total client satisfaction, the functions of this position are outlined below:

- Responsible to ensure that clients receive the proper sampling supplies.
- Accountable for response to client inquiries concerning sample status.
- Responsible for assistance to clients regarding the resolution of problems concerning COC.
- Ensuring that client specifications, when known, are met by communicating project and quality assurance requirements to the laboratory. Prepares a Quality Assurance Summary (QAS) or equivalent summary form as needed.
- Notifying the supervisors of incoming projects and sample delivery schedules.
- Accountable to clients for communicating sample progress in daily status meeting with agreed-upon due dates.
- Approves customer requested variances to methods and to standard laboratory protocols.
- Responsible for discussing with client any project-related problems, resolving service issues, and coordinating technical details with the laboratory staff.
- Reports client inquiries involving data quality issues or data acceptability to the facility, CSM, Laboratory Director, QA Manager, and to the appropriate staff.
- Responsible for staff familiarization with specific quotes, sample log-in review, and final report completeness.
- Monitor the status of all data package projects in-house to ensure timely and accurate delivery of reports. Reviews project data packages for completeness and compliance to client needs and have signature authority for final reports.
- Inform clients of data package-related problems and resolve service issues.
- Coordinate requests for sample containers and other services (data packages).
- Prepares re-issue requests for project data.
- Organizes bid activities for prospective new projects and clients.
- Has signature authority for laboratory reports.


#### 4.3 Deputies

The following table defines who assumes the responsibilities of key personnel in their absence:

Key Personnel	Deputy
Laboratory Director	QA Manager
QA Manager	Laboratory Director
Inorganic Technical Manager	Laboratory Director


Key Personnel	Deputy
Semivolatle Technical Manager	Laboratory Director
Volatile Technical Manager	Laboratory Director
EHS Coordinator	Laboratory Director
Project Manager Supervisor	Laboratory Director

Figure 4-1. Corporate and Laboratory Organization Charts

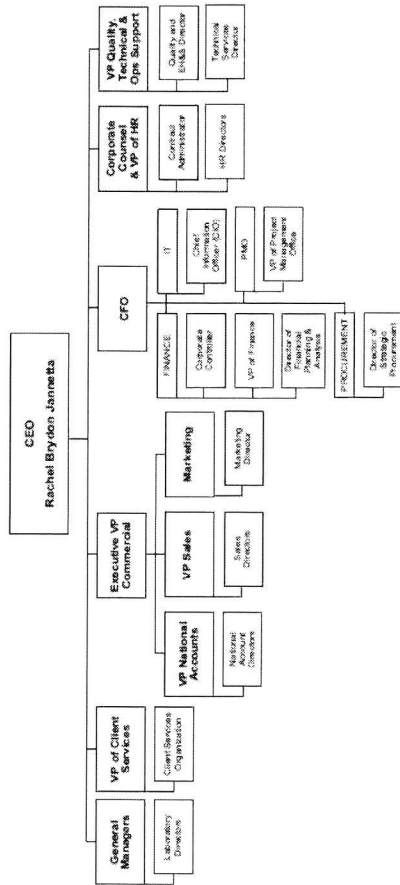

April 2014

### TestAmerica Seattle Organization Chart

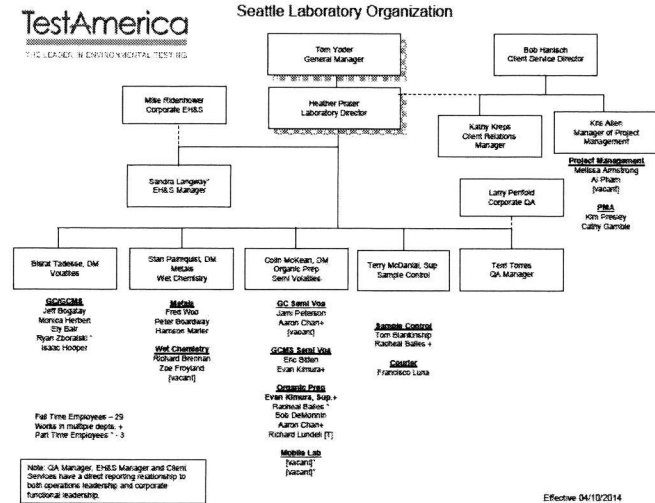

Effective 04/19/2014

## SECTION 5. QUALITY SYSTEM

### 5.1 Quality Policy Statement

It is TestAmerica's Policy to:

- ❖ Provide data of known quality to its clients by adhering to approved methodologies, regulatory requirements and the QA/QC protocols.
- ❖ Effectively manage all aspects of the laboratory and business operations by the highest ethical standards.
- ❖ Continually improve systems and provide support to quality improvement efforts in laboratory, administrative and managerial activities. TestAmerica recognizes that the implementation of a quality assurance program requires management's commitment and support as well as the involvement of the entire staff.
- ❖ Provide clients with the highest level of professionalism and the best service practices in the industry.
- ❖ To comply with the ISO/IEC 17025:2005(E) International Standard, the 2009 TNI Standard and to continually improve the effectiveness of the management system.

Every staff member at the laboratory plays an integral part in quality assurance and is held responsible and accountable for the quality of their work. It is, therefore, required that all laboratory personnel are trained and agree to comply with applicable procedures and requirements established by this document.

### 5.2 Ethics and Data Integrity

TestAmerica is committed to ensuring the integrity of its data and meeting the quality needs of its clients. The elements of TestAmerica's Ethics and Data Integrity Program include:

- An Ethics Policy (Corporate Policy No. CW-L-P-004) and Employee Ethics Statements.
- Ethics and Compliance Officers (ECOs).
- A Training Program.
- Self-governance through disciplinary action for violations.
- A Confidential mechanism for anonymously reporting alleged misconduct and a means for conducting internal investigations of all alleged misconduct. (Corporate SOP No. CW-L-S-002)
- Procedures and guidance for recalling data if necessary (Corporate SOP No. CW-L-S-002).
- Effective external and internal monitoring system that includes procedures for internal audits (Section 15).
- Produce results, which are accurate and include QA/QC information that meets client pre-defined Data Quality Objectives (DQOs).
- Present services in a confidential, honest and forthright manner.

- Provide employees with guidelines and an understanding of the Ethical and Quality Standards of our Industry.
- Operate our facilities in a manner that protects the environment and the health and safety of employees and the public.
- Obey all pertinent federal, state and local laws and regulations and encourage other members of our industry to do the same.
- Educate clients as to the extent and kinds of services available.
- Assert competency only for work for which adequate personnel and equipment are available and for which adequate preparation has been made.
- Promote the status of environmental laboratories, their employees, and the value of services rendered by them.

### 5.3 Quality System Documentation

The laboratory's Quality System is communicated through a variety of documents.

- Quality Assurance Manual – Each laboratory has a lab-specific quality assurance manual.
- Corporate SOPs and Policies – Corporate SOPs and Policies are developed for use by all relevant laboratories. They are incorporated into the laboratory's normal SOP distribution, training and tracking system. Corporate SOPs may be general or technical.
- Work Instructions – A subset of procedural steps, tasks or forms associated with an operation of a management system (e.g., checklists, preformatted bench sheets, forms).
- Laboratory SOPs – General and Technical
- Laboratory QA/QC Policy Memorandums

#### 5.3.1 Order of Precedence

In the event of a conflict or discrepancy between policies, the order of precedence is as follows:

- Corporate Quality Management Plan (CQMP)
- Corporate SOPs and Policies
- Laboratory QA/QC Policy Memorandum
- Laboratory Quality Assurance Manual (QAM)
- Laboratory SOPs and Policies
- Other (Work Instructions (WI), memos, flow charts, etc.)

Note: The laboratory has the responsibility and authority to operate in compliance with regulatory requirements of the jurisdiction in which the work is performed. Where the CQMP conflicts with those regulatory requirements, the regulatory requirements of the jurisdiction shall hold primacy. The laboratory's QAM shall take precedence over the CQMP in those cases.

#### 5.4 QA/QC Objectives for the Measurement of Data

Quality Assurance (QA) and Quality Control (QC) are activities undertaken to achieve the goal of producing data that accurately characterize the sites or materials that have been sampled. Quality Assurance is generally understood to be more comprehensive than Quality Control. Quality Assurance can be defined as the integrated system of activities that ensures that a product or service meets defined standards.

Quality Control is generally understood to be limited to the analyses of samples and to be synonymous with the term "analytical quality control". QC refers to the routine application of statistically based procedures to evaluate and control the accuracy of results from analytical measurements. The QC program includes procedures for estimating and controlling precision and bias and for determining reporting limits.

Request for Proposals (RFPs) and Quality Assurance Project Plans (QAPP) provide a mechanism for the client and the laboratory to discuss the data quality objectives in order to ensure that analytical services closely correspond to client needs. The client is responsible for developing the QAPP. In order to ensure the ability of the laboratory to meet the Data Quality Objectives (DQOs) specified in the QAPP, clients are advised to allow time for the laboratory to review the QAPP before being finalized. Additionally, the laboratory will provide support to the client for developing the sections of the QAPP that concern laboratory activities.

Historically, laboratories have described their QC objectives in terms of precision, accuracy, representativeness, comparability, completeness, selectivity and sensitivity (PARCCSS).

##### 5.4.1 Precision

The laboratory objective for precision is to meet the performance for precision demonstrated for the methods on similar samples and to meet data quality objectives of the EPA and/or other regulatory programs. Precision is defined as the degree of reproducibility of measurements under a given set of analytical conditions (exclusive of field sampling variability). Precision is documented on the basis of replicate analysis, usually duplicate or matrix spike (MS) duplicate samples.

##### 5.4.2 Accuracy

The laboratory objective for accuracy is to meet the performance for accuracy demonstrated for the methods on similar samples and to meet data quality objectives of the EPA and/or other regulatory programs. Accuracy is defined as the degree of bias in a measurement system. Accuracy may be documented through the use of laboratory control samples (LCS) and/or MS. A statement of accuracy is expressed as an interval of acceptance recovery about the mean recovery.

##### 5.4.3 Representativeness

The laboratory objective for representativeness is to provide data which is representative of the sampled medium. Representativeness is defined as the degree to which data represent a characteristic of a population or set of samples and is a measurement of both analytical and field sampling precision. The representativeness of the analytical data is a function of the

procedures used in procuring and processing the samples. The representativeness can be documented by the relative percent difference between separately procured, but otherwise identical samples or sample aliquots.

The representativeness of the data from the sampling sites depends on both the sampling procedures and the analytical procedures. The laboratory may provide guidance to the client regarding proper sampling and handling methods in order to assure the integrity of the samples.

##### 5.4.4 Comparability

The comparability objective is to provide analytical data for which the accuracy, precision, representativeness and reporting limit statistics are similar to these quality indicators generated by other laboratories for similar samples, and data generated by the laboratory over time.

The comparability objective is documented by inter-laboratory studies carried out by regulatory agencies or carried out for specific projects or contracts, by comparison of periodically generated statements of accuracy, precision and reporting limits with those of other laboratories.

##### 5.4.5 Completeness

The completeness objective for data is 90% (or as specified by a particular project), expressed as the ratio of the valid data to the total data over the course of the project. Data will be considered valid if they are adequate for their intended use. Data usability will be defined in a QAPP, project scope or regulatory requirement. Data validation is the process for reviewing data to determine its usability and completeness. If the completeness objective is not met, actions will be taken internally and with the data user to improve performance. This may take the form of an audit to evaluate the methodology and procedures as possible sources for the difficulty or may result in a recommendation to use a different method.

##### 5.4.6 Selectivity

Selectivity is defined as: The capability of a test method or instrument to respond to a target substance or constituent in the presence of non-target substances. Target analytes are separated from non-target constituents and subsequently identified/detected through one or more of the following, depending on the analytical method: extractions (separation), digestions (separation), interelement corrections (separation), use of matrix modifiers (separation), specific retention times (separation and identification), confirmations with different columns or detectors (separation and identification), specific wavelengths (identification), specific mass spectra (identification), specific electrodes (separation and identification), etc..

##### 5.4.7 Sensitivity

Sensitivity refers to the amount of analyte necessary to produce a detector response that can be reliably detected (Method Detection Limit) or quantified (Reporting Limit).

### 5.5 Criteria for Quality Indicators

The laboratory maintains quality control limits in the Laboratory Information Management System that summarize the precision and accuracy acceptability limits for performed analyses. This summary includes an effective date, is updated each time new limits are generated and are managed by the laboratory's QA department. Unless otherwise noted, limits within these tables are laboratory generated. Some acceptability limits are derived from US EPA methods when they are required. Where US EPA method limits are not required, the laboratory has developed limits from evaluation of data from similar matrices. Criteria for development of control limits is contained in SOP TA-QA-0600 Control Charting and Establishing Method Warning and Action Limits.

### 5.6 Statistical Quality Control

Statistically-derived precision and accuracy limits are required by selected methods (such as SW-846) and programs [such as the Ohio Voluntary Action Plan (VAP)]. The laboratory routinely utilizes statistically-derived limits to evaluate method performance and determine when corrective action is appropriate. The analysts are instructed to use the current limits in the laboratory (dated and approved by the Technical Manager and QA Manager) and entered into the Laboratory Information Management System (LIMS). The Quality Assurance department maintains an archive of all limits used within the laboratory via the LIMS. If a method defines the QC limits, the method limits are used.

If a method requires the generation of historical limits, the lab develops such limits from recent data in the QC database of the LIMS following the guidelines described in Section 24. All calculations and limits are documented and dated when approved and effective. On occasion, a client requests contract-specified limits for a specific project.

Current QC limits are entered and maintained in the LIMS analyte database. As sample results and the related QC are entered into LIMS, the sample QC values are compared with the limits in LIMS to determine if they are within the acceptable range. The analyst then evaluates if the sample needs to be rerun or re-extracted/rerun or if a comment should be added to the report explaining the reason for the QC outlier.

#### 5.6.1 QC Charts

As the QC limits are calculated, QC charts are generated showing warning and control limits for the purpose of evaluating trends. The QA Manager evaluates these to determine if adjustments need to be made or for corrective actions to methods. All findings are documented and kept on file. (Refer to SOP TA-QA-0600 Control Charting and Establishing Method Warning and Action Limits.)

### 5.7 Quality System Metrics

In addition to the QC parameters discussed above, the entire Quality System is evaluated on a monthly basis through the use of specific metrics (refer to Section 16). These metrics are used to drive continuous improvement in the laboratory's Quality System.

## SECTION 6. DOCUMENT CONTROL

### 6.1 Overview

The QA Department is responsible for the control of documents used in the laboratory to ensure that approved, up-to-date documents are in circulation and out-of-date (obsolete) documents are archived or destroyed. The following documents, at a minimum, must be controlled:

- Laboratory Quality Assurance Manual
- Laboratory Standard Operating Procedures (SOP)
- Laboratory Policies
- Work Instructions and Forms
- Corporate Policies and Procedures distributed outside the intranet

Corporate Quality posts Corporate Manuals, SOPs, Policies, Work Instructions, White Papers and Training Materials on the company intranet site. These Corporate documents are only considered controlled when they are read on the intranet site. Printed copies are considered uncontrolled unless the laboratory physically distributes them as controlled documents. A detailed description of the procedure for issuing, authorizing, controlling, distributing, and archiving Corporate documents is found in Corporate SOP No. CW-Q-S-001, Corporate Document Control and Archiving. The laboratory's internal document control procedures are defined in SOPs TA-QA-0528, (Document Control) and TA-QA-0506 (Document Archiving).

The laboratory QA Department also maintains access to various references and document sources integral to the operation of the laboratory. This includes reference methods and regulations. Instrument manuals (hard or electronic copies) are also maintained by the laboratory.

The laboratory maintains control of records for raw analytical data and supporting records such as audit reports and responses, logbooks, standard logs, training files, MDL studies, Proficiency Testing (PT) studies, certifications and related correspondence, and nonconformance memos/corrective action reports. Raw analytical data consists of bound logbooks, instrument printouts, any other notes, magnetic media, electronic data and final reports.

### 6.2 Document Approval and Issue

The pertinent elements of a document control system for each document include a unique document title and number, pagination, the total number of pages of the item or an 'end of document' page, the effective date, revision number and the laboratory's name. The QA personnel are responsible for the maintenance of this system.

Controlled documents are authorized by the QA Department. In order to develop a new document, a technical manager submits an electronic draft to the QA Department for suggestions and approval before use. Upon approval, QA personnel add the identifying version information to the document and retains that document as the official document on file. That document is then provided to all applicable operational units (may include electronic access). Controlled documents are identified as such and records of their distribution are kept by the QA Department. Document control may be achieved by either electronic or hardcopy distribution.



## Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 7  
Effective Date: 4/23/2014  
Page 33 of 165

The QA Department maintains a list of the official versions of controlled documents.

Quality System Policies and Procedures will be reviewed at a minimum of annually and revised as appropriate. Changes to documents occur when a procedural change warrants.

### 6.3 Procedures for Document Control Policy

For changes to the QA Manual, refer to SOP No. TA-QA-0528 (Document Control). Uncontrolled copies must not be used within the laboratory. Previous revisions and back-up data are stored by the QA department. Electronic copies are stored on the Public server in the QAM folder for the applicable revision.

For changes to SOPs, refer to SOP No. TA-QA-0500, Standard Operating Procedures. The SOP identified above also defines the process of changes to SOPs.

Forms, worksheets, work instructions and information are organized by department in the QA office. There is a table of contents. Electronic versions are kept on a hard drive in the QA department; hard copies are kept in QA files. The procedure for the care of these documents is in SOP TA-QA-0528 (Document Control).

### 6.4 Obsolete Documents

All invalid or obsolete documents are removed, or otherwise prevented from unintended use. The laboratory has specific procedures as described above to accomplish this. In general, obsolete documents are collected from employees according to distribution lists and are marked obsolete on the cover or destroyed. At least one copy of the obsolete document is archived according to SOP No. TA-QA-0528 (Document Control).

## SECTION 7. SERVICE TO THE CLIENT

### 7.1 Overview

The laboratory has established procedures for the review of work requests and contracts, oral or written. The procedures include evaluation of the laboratory's capability and resources to meet the contract's requirements within the requested time period. All requirements, including the methods to be used, must be adequately defined, documented and understood. For many environmental sampling and analysis programs, testing design is site or program specific and does not necessarily "fit" into a standard laboratory service or product. It is the laboratory's intent to provide both standard and customized environmental laboratory services to our clients.

A thorough review of technical and QC requirements contained in contracts is performed to ensure project success. The appropriateness of requested methods, and the lab's capability to perform them must be established. Projects, proposals and contracts are reviewed for adequately defined requirements and the laboratory's capability to meet those requirements. Alternate test methods that are capable of meeting the clients' requirements may be proposed by the lab. A review of the lab's capability to analyze non-routine analytes is also part of this review process.

## Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 7  
Effective Date: 4/23/2014  
Page 34 of 165

All projects, proposals and contracts are reviewed for the client's requirements in terms of compound lists, test methodology requested, sensitivity (detection and reporting levels), accuracy, and precision requirements (% Recovery and RPD). The reviewer ensures that the laboratory's test methods are suitable to achieve these requirements and that the laboratory holds the appropriate certifications and approvals to perform the work. The laboratory and any potential subcontract laboratories must be certified, as required, for all proposed tests.

The laboratory must determine if it has the necessary physical, personnel and information resources to meet the contract, and if the personnel have the expertise needed to perform the testing requested. Each proposal is checked for its impact on the capacity of the laboratory's equipment and personnel. As part of the review, the proposed turnaround time will be checked for feasibility.

Electronic or hard copy deliverable requirements are evaluated against the laboratory's capacity for production of the documentation.

If the laboratory cannot provide all services but intends to subcontract such services, whether to another TestAmerica facility or to an outside firm, this will be documented and discussed with the client prior to contract approval. (Refer to Section 8 for Subcontracting Procedures.)

The laboratory informs the client of the results of the review if it indicates any potential conflict, deficiency, lack of accreditation, or inability of the lab to complete the work satisfactorily. Any discrepancy between the client's requirements and the laboratory's capability to meet those requirements is resolved in writing before acceptance of the contract. It is necessary that the contract be acceptable to both the laboratory and the client. Amendments initiated by the client and/or TestAmerica, are documented in writing.

All contracts, QAPPs, Sampling and Analysis Plans (SAPs), contract amendments, and documented communications become part of the project record.

The same contract review process used for the initial review is repeated when there are amendments to the original contract by the client, and the participating personnel are informed of the changes.

### 7.2 Review Sequence and Key Personnel

Appropriate personnel will review the work request at each stage of evaluation.

For routine projects and other simple tasks, a review by the Project Manager (PM) is considered adequate. The PM confirms that the laboratory has any required certifications, that it can meet the clients' data quality and reporting requirements and that the lab has the capacity to meet the clients turn around needs. It is recommended that, where there is a sales person assigned to the account, an attempt should be made to contact that sales person to inform them of the incoming samples.

For new, complex or large projects, the proposed contract is given to the Sales Directors, who will decide which lab will receive the work based on the scope of work and other requirements, including certification, testing methodology, and available capacity to perform the work. The

contract review process is outlined in TestAmerica's Corporate SOP No. CA-L-P-002, Contract Compliance Policy.

This review encompasses all facets of the operation. The scope of work is distributed to the appropriate personnel, as needed based on scope of contract, to evaluate all of the requirements shown above (not necessarily in the order below):

- Legal & Contracts Director
- General Manager
- Customer Services Manager
- Laboratory Operations Manager
- Laboratory and/or Corporate Technical Managers / Directors
- Corporate Information Technology Managers
- Account Executives
- Laboratory and/or Corporate Quality
- Laboratory Environmental Health and Safety Managers
- The Laboratory Director reviews the formal laboratory quote and makes final acceptance for their facility.

The Sales Director, Legal Contracts Director or Account Executive then submits the final proposal to the client.

In the event that one of the above personnel is not available to review the contract, his or her back-up will fulfill the review requirements.

The Legal & Contracts Director maintains copies of all signed contracts. The Laboratory Director also maintains copies of all signed contracts on site.

### 7.3 Documentation

Appropriate records are maintained for every contract or work request. All stages of the contract review process are documented and include records of any significant changes. The reviewed contracts or QAPPs, along with records identifying issues or changes (at a minimum the Additional Project Notes from the quote) will be scanned and stored on the Public server in the QA folder.

The contract will be distributed to and maintained by the appropriate sales/marketing personnel and the Account Executive. A copy of the contract and formal quote will be filed with the laboratory PM and the Laboratory Director.

Records are maintained of pertinent discussions with a client relating to the client's requirements or the results of the work during the period of execution of the contract. The PM keeps a phone log of conversations with the client and all actions agreed upon between the PM and the client are officially documented with a follow-up email.

### 7.3.1 Project-Specific Quality Planning

Communication of contract specific technical and QC criteria is an essential activity in ensuring the success of site specific testing programs. To achieve this goal, the laboratory assigns a PM to each client. It is the PM's responsibility to ensure that project-specific technical and QC requirements are effectively evaluated and communicated to the laboratory personnel before and during the project. QA department involvement may be needed to assist in the evaluation of custom QC requirements.

PM's are the primary client contact and they ensure resources are available to meet project requirements. Although PM's do not have direct reports or staff in production, they coordinate opportunities and work with laboratory management and supervisory staff to ensure available resources are sufficient to perform work for the client's project. Project management is positioned between the client and laboratory resources.

Prior to work on a new project, the dissemination of project information and/or project opening meetings may occur to discuss schedules and unique aspects of the project. Items to be discussed may include the project technical profile, turnaround times, holding times, methods, analyte lists, reporting limits, deliverables, sample hazards, or other special requirements. The PM introduces new projects to the laboratory staff through project kick-off meetings or to the supervisory staff during production meetings. These meetings provide direction to the laboratory staff in order to maximize production and client satisfaction, while maintaining quality. In addition, project notes may be associated with each sample batch as a reminder upon sample receipt and analytical processing.

During the project, any change that may occur within an active project is agreed upon between the client/regulatory agency and the PM/laboratory. These changes (e.g., use of a non-standard method or modification of a method) and approvals must be documented prior to implementation. Documentation pertains to any document, e.g., letter, e-mail, variance, contract addendum, which has been signed by both parties.

Such changes are communicated to the laboratory during production or status meetings. These changes are also updated in the project notes and are introduced to the managers at these meetings. The laboratory staff is then introduced to the modified requirements via the PM or the individual laboratory Technical Manager. After the modification is implemented into the laboratory process, documentation of the modification is made in the case narrative of the data report(s).

The laboratory strongly encourages client visits to the laboratory and for formal/informal information sharing session with employees in order to effectively communicate ongoing client needs as well as project specific details for customized testing programs.

### 7.4 Special Services

The laboratory cooperates with clients and their representatives to monitor the laboratory's performance in relation to work performed for the client. It is the laboratory's goal to meet all client requirements in addition to statutory and regulatory requirements. The laboratory has procedures to ensure confidentiality to clients (Section 15 and 25).

**Note:** ISO/IEC 17025 states that a laboratory "shall afford clients or their representatives cooperation to clarify the client's request". This topic is discussed in Section 7.

The laboratory's standard procedures for reporting data are described in Section 25. Special services are also available and provided upon request. These services include:

- Reasonable access for our clients or their representatives to the relevant areas of the laboratory for the witnessing of tests performed for the client.
- Assist client-specified third party data validators as specified in the client's contract.
- Supplemental information pertaining to the analysis of their samples. **Note:** An additional charge may apply for additional data/information that was not requested prior to the time of sample analysis or previously agreed upon.

#### 7.5 Client Communication

Project managers are the primary communication link to the clients. They shall inform their clients of any delays in project completion as well as any non-conformances in either sample receipt or sample analysis. Project management will maintain ongoing client communication throughout the entire client project.

Technical Managers are available to discuss any technical questions or concerns that the client may have.

#### 7.6 Reporting

The laboratory works with our clients to produce any special communication reports required by the contract.

#### 7.7 Client Surveys

The laboratory assesses both positive and negative client feedback. The results are used to improve overall laboratory quality and client service. TestAmerica's Sales and Marketing teams periodically develop lab- and client-specific surveys to assess client satisfaction.

### SECTION 8. SUBCONTRACTING OF TESTS

#### 8.1 Overview

For the purpose of this quality manual, the phrase subcontract laboratory refers to a laboratory external to the TestAmerica laboratories. The phrase "work sharing" refers to internal transfers of samples between the TestAmerica laboratories. The term outsourcing refers to the act of subcontracting tests.

When contracting with our clients, the laboratory makes commitments regarding the services to be performed and the data quality for the results to be generated. When the need arises to outsource testing for our clients because project scope, changes in laboratory capabilities, capacity or unforeseen circumstances, we must be assured that the subcontractors or work sharing laboratories understand the requirements and will meet the same commitments we

have made to the client. Refer to TestAmerica's Corporate SOP's on Subcontracting Procedures (CA-L-S-002) and the Work Sharing Process (CA-C-S-001).

When outsourcing analytical services, the laboratory will assure, to the extent necessary, that the subcontract or work sharing laboratory maintains a program consistent with the requirements of this document, the requirements specified in TNI/ISO 17025 and/or the client's Quality Assurance Project Plan (QAPP). All QC guidelines specific to the client's analytical program are transmitted to the subcontractor and agreed upon before sending the samples to the subcontract facility. Additionally, work requiring accreditation will be placed with an appropriately accredited laboratory. The laboratory performing the subcontracted work will be identified in the final report, as will non-NELAC accredited work where required.

Project Managers (PMs), Project Manager Supervisor, or Account Executives (AE) for the Export Lab are responsible for obtaining client approval prior to outsourcing any samples. The laboratory will advise the client of a subcontract or work sharing arrangement in writing and when possible approval from the client shall be retained in the project folder.

**Note:** In addition to the client, some regulating agencies (e.g., USDA) or contracts (e.g., certain USACE projects) may require notification prior to placing such work.

#### 8.2 Qualifying and Monitoring Subcontractors

Whenever a PM or Customer Service Manager becomes aware of a client requirement or laboratory need where samples must be outsourced to another laboratory, the other laboratory(s) shall be selected based on the following:

- The first priority is to attempt to place the work in a qualified TestAmerica laboratory;
- Firms specified by the client for the task (Documentation that a subcontractor was designated by the client must be maintained with the project file. This documentation can be as simple as placing a copy of an e-mail from the client in the project folder);
- Firms listed as pre-qualified and currently under a subcontract with TestAmerica: A listing of all approved subcontracting laboratories is available on the TestAmerica intranet site. Supporting documentation is maintained by corporate offices and by the TestAmerica laboratory originally requesting approval of the subcontract lab. Verify necessary accreditation, where applicable, (e.g., on the subcontractors TNI, A2LA accreditation or State Certification).
- Firms identified in accordance with the company's Small Business Subcontracting program as small, women-owned, veteran-owned and/or minority-owned businesses;
- NELAC or A2LA accredited laboratories.
- In addition, the firm must hold the appropriate certification to perform the work required.

All TestAmerica laboratories are pre-qualified for work sharing provided they hold the appropriate accreditations, can adhere to the project/program requirements, and the client approved sending samples to that laboratory. The client must provide acknowledgement that the samples can be sent to that facility (an e-mail is sufficient documentation or if acknowledgement is verbal, the date, time, and name of person providing acknowledgement must be documented). The originating laboratory is responsible for communicating all technical,

quality, and deliverable requirements as well as other contract needs. (Corporate SOP No. CA-C-S-001, Work Sharing Process).

When the potential sub-contract laboratory has not been previously approved, PMs may nominate a laboratory as a subcontractor based on need. The decision to nominate a laboratory must be approved by the Laboratory Director. The Laboratory Director requests that the QA Manager begin the process of approving the subcontract laboratory as outlined in Corporate SOP No. CA-L-S-002, Subcontracting Procedures. The client must provide acknowledgement that the samples can be sent to that facility (an e-mail is sufficient documentation or if acknowledgement is verbal, the date, time, and name of person providing acknowledgement must be documented).

**8.2.1** Once the appropriate accreditation and legal information is received by the laboratory, it is evaluated for acceptability (where applicable) and forwarded to Corporate Contracts for formal contracting with the laboratory. They will add the lab to the approved list on the intranet site and notify the finance group for JD Edwards.

**8.2.2** The client will assume responsibility for the quality of the data generated from the use of a subcontractor they have requested the lab to use. The qualified subcontractors on the intranet site are known to meet minimal standards. TestAmerica does not certify laboratories. The subcontractor is on our approved list and can only be recommended to the extent that we would use them.

**8.2.3** The status and performance of qualified subcontractors will be monitored periodically by the Corporate Contracts and/or Quality Departments. Any problems identified will be brought to the attention of TestAmerica's Corporate Finance or Corporate Quality personnel.

- Complaints shall be investigated. Documentation of the complaint, investigation and corrective action will be maintained in the subcontractor's file on the intranet site. Complaints are posted using the Vendor Performance Report.
- Information shall be updated on the intranet when new information is received from the subcontracted laboratories.
- Subcontractors in good standing will be retained on the intranet listing. The QA Manager will notify all TestAmerica laboratories, Corporate Quality and Corporate Contracts if any laboratory requires removal from the intranet site. This notification will be posted on the intranet site and e-mailed to all Laboratory Directors, QA Managers and Sales Personnel.

### **8.3 Oversight and Reporting**

The PM must request that the selected subcontractor be presented with a subcontract, if one is not already executed between the laboratory and the subcontractor. The subcontract must include terms which flow down the requirements of our clients, either in the subcontract itself or through the mechanism of work orders relating to individual projects. A standard subcontract and the Lab Subcontractor Vendor Package (posted on the intranet) can be used to accomplish this, and the Legal & Contracts Director can tailor the document or assist with negotiations, if needed. The PM or CSM responsible for the project must advise and obtain client consent to the subcontract as appropriate, and provide the scope of work to ensure that the proper requirements are made a part of the subcontract and are made known to the subcontractor.

Prior to sending samples to the subcontracted laboratory, the PM confirms their certification status to determine if it's current and scope-inclusive. The information is documented on a Subcontracted Sample Form (Figure 8-1) and the form is retained in the project folder. For TestAmerica laboratories, certifications can be viewed on the company's TotalAccess Database.

The Sample Control department is responsible for ensuring compliance with QA requirements and applicable shipping regulations when shipping samples to a subcontracted laboratory.

All subcontracted samples must be accompanied by a TestAmerica Chain of Custody (COC). A copy of the original COC sent by the client must also be included with all samples workshared within TestAmerica. Client CoCs are only forwarded to external subcontractors when samples are shipped directly from the project site to the subcontractor lab. Under routine circumstances, client CoCs are not provided to external subcontractors.

Through communication with the subcontracted laboratory, the PM monitors the status of the subcontracted analyses, facilitates successful execution of the work, and ensures the timeliness and completeness of the analytical report.

Non-NELAC accredited work must be identified in the subcontractor's report as appropriate. If NELAC accreditation is not required, the report does not need to include this information.

Reports submitted from subcontractor laboratories are not altered and are included in their original form in the final project report. This clearly identifies the data as being produced by a subcontractor facility. If subcontract laboratory data is incorporated into the laboratories EDD (i.e., imported), the report must explicitly indicate which lab produced the data for which methods and samples.

**Note:** The results submitted by a TestAmerica work sharing laboratory may be transferred electronically and the results reported by the TestAmerica work sharing lab are identified on the final report. The report must explicitly indicate which lab produced the data for which methods and samples. The final report must include a copy of the completed COC for all work sharing reports.

### **8.4 Contingency Planning**

The Laboratory Director may waive the full qualification of a subcontractor process temporarily to meet emergency needs. However, this decision & justification must be documented in the project files, and the 'Purchase Order Terms and Conditions For Subcontracted Laboratory Services' must be sent with the samples and Chain-of-Custody. In the event this provision is utilized, the laboratory (e.g., PM) will be required to verify and document the applicable accreditations of the subcontractor. All other quality and accreditation requirements will still be applicable, but the subcontractor need not have signed a subcontract with TestAmerica at this time. The comprehensive approval process must then be initiated within 30 calendar days of subcontracting.

Figure 8-1.

**Example - Subcontracted Sample Form**

**Date/Time:** \_\_\_\_\_

**Subcontracted Laboratory Information:**

- Subcontractor's Name: \_\_\_\_\_
- Subcontractor Point of Contact: \_\_\_\_\_
- Subcontractor's Address: \_\_\_\_\_
- Subcontractor's Phone: \_\_\_\_\_
- Analyte/Method: \_\_\_\_\_
- Certified for State of Origin: \_\_\_\_\_
- NELAC Certified: \_\_\_\_\_  
Yes \_\_\_\_\_ No \_\_\_\_\_
- USDA Permit ( \_\_ Domestic \_\_ Foreign) Yes \_\_\_\_\_ No \_\_\_\_\_
- A2LA (or ISO 17025) Certified: Yes \_\_\_\_\_ No \_\_\_\_\_
- CLP-like Required: Yes \_\_\_\_\_ No \_\_\_\_\_  
(Full doc required)
- Requested Sample Due Date: \_\_\_\_\_  
(Must be put on COC)

**Project Manager:** \_\_\_\_\_

**Laboratory Sample # Range:** \_\_\_\_\_  
(Only of Subcontracted Samples)

**Laboratory Project Number (Billing Control #):** \_\_\_\_\_

All subcontracted samples are to be sent via bonded carrier and Priority Overnight. Please attach tracking number below and maintain these records in the project files.

**PM Signature** \_\_\_\_\_ **Date** \_\_\_\_\_

**SECTION 9. PURCHASING SERVICES AND SUPPLIES****9.1 Overview**

Evaluation and selection of suppliers and vendors is performed, in part, on the basis of the quality of their products, their ability to meet the demand for their products on a continuous and short term basis, the overall quality of their services, their past history, and competitive pricing. This is achieved through evaluation of objective evidence of quality furnished by the supplier, which can include certificates of analysis, recommendations, and proof of historical compliance with similar programs for other clients. To ensure that quality critical consumables and equipment conform to specified requirements, which may affect quality, all purchases from specific vendors are approved by a member of the supervisory or management staff. Capital expenditures are made in accordance with TestAmerica's Corporate Controlled Purchases Procedure, SOP No. CW-F-S-007.

Contracts will be signed in accordance with TestAmerica's Corporate Authorization Matrix Policy, Policy No. CW-F-P-002. Request for Proposals (RFP's) will be issued where more information is required from the potential vendors than just price. Process details are available in TestAmerica's Corporate Procurement and Contracts Policy (Policy No. CW-F-P-004). RFP's allow TestAmerica to determine if a vendor is capable of meeting requirements such as supplying all of the TestAmerica facilities, meeting required quality standards and adhering to necessary ethical and environmental standards. The RFP process also allows potential vendors to outline any additional capabilities they may offer.

**9.2 Glassware**

Glassware used for volumetric measurements must be Class A or verified for accuracy according to laboratory procedure. Pyrex (or equivalent) glass should be used where possible. For safety purposes, thick-wall glassware should be used where available.

**9.3 Reagents, Standards & Supplies**

Purchasing guidelines for equipment and reagents must meet the requirements of the specific method and testing procedures for which they are being purchased. Solvents and acids are pre-tested in accordance with TestAmerica's Corporate SOP on Solvent & Acid Lot Testing & Approval, SOP No. CA-Q-S-001.

**9.3.1 Purchasing**

Chemical reagents, solvents, glassware, and general supplies are ordered as needed to maintain sufficient quantities on hand. Materials used in the analytical process must be of a known quality. The wide variety of materials and reagents available makes it advisable to specify recommendations for the name, brand, and grade of materials to be used in any determination. This information is contained in the method SOP. The analyst completes the Material Request Sheet when requesting reagents, standards, or supplies that are not stocked in on-site consignment system that contains items approved for laboratory use.

The analyst must provide the master item number (from the master item list that has been approved by the Technical Manager), item description, package size, catalogue page number,



## Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 7  
Effective Date: 4/23/2014  
Page 43 of 165

and the quantity needed. If an item being ordered is not the exact item requested, approval must be obtained from the Technical Manager prior to placing the order. The authorized requisitioning agent places the order.

### 9.3.2 Receiving

It is the responsibility of the receiving clerk to receive the shipment. It is the responsibility of the analyst who ordered the materials to document the date materials were received. Once the ordered reagents or materials are received, the analyst compares the information on the label or packaging to the original order to ensure that the purchase meets the quality level specified. Material Safety Data Sheets (MSDSs) are available online through the Company's intranet website. Anyone may review these for relevant information on the safe handling and emergency precautions of on-site chemicals.

### 9.3.3 Specifications

Methods in use in the laboratory specify the grade of reagent that must be used in the procedure. If the quality of the reagent is not specified, analytical reagent grade will be used. It is the responsibility of the analyst to check the procedure carefully for the suitability of grade of reagent.

Chemicals must not be used past the manufacturer's expiration date and must not be used past the expiration time noted in a method SOP. If expiration dates are not provided, the laboratory may contact the manufacturer to determine an expiration date.

The laboratory assumes a five year expiration date on inorganic dry chemicals and solvents unless noted otherwise by the manufacturer or by the reference source method. Chemicals/solvents should not be used past the manufacturer's or SOPs expiration date unless 'verified' (refer to item 3 listed below).

- An expiration date **cannot** be extended if the dry chemical/solvent is discolored or appears otherwise physically degraded, the dry chemical/solvent must be discarded.
- Expiration dates can be extended if the dry chemical/solvent is found to be satisfactory based on acceptable performance of quality control samples (Continuing Calibration Verification (CCV), Blanks, Laboratory Control Sample (LCS), etc.).
- If the dry chemical/solvent is used for the preparation of standards, the expiration dates can be extended 6 months if the dry chemical/solvent is compared to an unexpired independent source in performing the method and the performance of the dry chemical/solvent is found to be satisfactory. The comparison must show that the dry chemical/solvent meets CCV limits. The comparison studies are maintained in the QA office.

Wherever possible, standards must be traceable to national or international standards of measurement or to national or international reference materials. Records to that effect are available to the user.

Compressed gases in use are checked for pressure and secure positioning daily. To prevent a tank from going to dryness or introducing potential impurities, the pressure should be closely watched as it decreases to approximately 15% of the original reading, at which point it should

## Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 7  
Effective Date: 4/23/2014  
Page 44 of 165

be replaced. For example, a standard sized laboratory gas cylinder containing 3,000 psig of gas should be replaced when it drops to approximately 500 psig. The quality of the gases must meet method or manufacturer specification or be of a grade that does not cause any analytical interference.

Water used in the preparation of standards or reagents must have a specific conductivity of less than 1-  $\mu\text{mho/cm}$  (or specific resistivity of greater than 1.0 megohm-cm) at 25°C. The specific conductivity is checked and recorded daily. If the water's specific conductivity is greater than the specified limit, the Facility Manager and appropriate Technical Manager must be notified immediately in order to notify all departments, decide on cessation (based on intended use) of activities, and make arrangements for correction.

The laboratory may purchase reagent grade (or other similar quality) water for use in the laboratory. This water must be certified "clean" by the supplier for all target analytes or otherwise verified by the laboratory prior to use. This verification is documented.

Standard lots are verified before first time use if the laboratory switches manufacturers or has historically had a problem with the type of standard.

Purchased bottleware used for sampling must be certified clean and the certificates must be maintained. If uncertified sampling bottleware is purchased, all lots must be verified clean prior to use. This verification must be maintained.

Records of manufacturer's certification and traceability statements are stored as scanned images in the LIMS Reagent module. These records include date of receipt, lot number (when applicable), and expiration date (when applicable). Incorporation of the item into the record indicates that the analyst has compared the new certificate with the previous one for the same purpose and that no difference is noted, unless approved and so documented by the Technical Manager or QA Manager.

### 9.3.4 Storage

Reagent and chemical storage is important from the aspects of both integrity and safety. Light-sensitive reagents may be stored in brown-glass containers. Storage conditions are per the Corporate Environmental Health & Safety Manual (Corp. Doc. No. CW-E-M-001) and method SOPs or manufacturer instructions.

### 9.4 Purchase of Equipment / Instruments / Software

When a new piece of equipment is needed, either for additional capacity or for replacing inoperable equipment, the analyst or supervisor makes a supply request to the Technical Manager and/or the Laboratory Director. If they agree with the request, the procedures outlined in TestAmerica's Corporate Policy No. CA-T-P-001, Qualified Products List, are followed. A decision is made as to which piece of equipment can best satisfy the requirements. The appropriate written requests are completed and purchasing places the order.

Upon receipt of a new or used piece of equipment, an identification name is assigned and added to the equipment list. IT must also be notified so that they can synchronize the instrument for back-ups. Its capability is assessed to determine if it is adequate or not for the

specific application. For instruments, a calibration curve is generated, followed by MDLs, Demonstration of Capabilities (DOCs), and other relevant criteria (refer to Section 19). For software, its operation must be deemed reliable and evidence of instrument verification must be retained by the QA Department. Software certificates supplied by the vendors are filed with the LIMS Administrator. The manufacturer's operation manual is retained at the bench.

#### 9.5 Services

Service to analytical instruments (except analytical balances) is performed on an as needed basis. Routine preventative maintenance is discussed in Section 20. The need for service is determined by analysts and/or Technical Managers. The service providers that perform the services are approved by the Technical Manager.

#### 9.6 Suppliers

TestAmerica selects vendors through a competitive proposal / bid process, strategic business alliances or negotiated vendor partnerships (contracts). This process is defined in the Corporate Finance documents on Vendor Selection (SOP No. CW-F-S-018) and Procurement & Contracts Policy (Policy No. CW-F-P-004). The level of control used in the selection process is dependent on the anticipated spending amount and the potential impact on TestAmerica business. Vendors that provide test and measuring equipment, solvents, standards, certified containers, instrument related service contracts or subcontract laboratory services shall be subject to more rigorous controls than vendors that provide off-the-shelf items of defined quality that meet the end use requirements. The JD Edwards purchasing system includes all suppliers/vendors that have been approved for use.

Evaluation of suppliers is accomplished by ensuring the supplier ships the product or material ordered and that the material is of the appropriate quality. This is documented by signing off on packing slips or other supply receipt documents. The purchasing documents contain the data that adequately describe the services and supplies ordered.

Any issues of vendor performance are to be reported immediately by the laboratory staff to the Corporate Purchasing Group by completing a Vendor Performance Report.

The Corporate Purchasing Group will work through the appropriate channels to gather the information required to clearly identify the problem and will contact the vendor to report the problem and to make any necessary arrangements for exchange, return authorization, credit, etc.

As deemed appropriate, the Vendor Performance Reports will be summarized and reviewed to determine corrective action necessary, or service improvements required by vendors.

The laboratory has access to a listing of all approved suppliers of critical consumables, supplies and services. This information is provided through the JD Edwards purchasing system.

#### 9.6.1 New Vendor Procedure

TestAmerica employees who wish to request the addition of a new vendor must complete a J.D. Edwards Vendor Add Request Form.

New vendors are evaluated based upon criteria appropriate to the products or services provided as well as their ability to provide those products and services at a competitive cost. Vendors are also evaluated to determine if there are ethical reasons or potential conflicts of interest with TestAmerica employees that would make it prohibitive to do business with them as well as their financial stability. The QA Department and/or the Technology Manager are consulted with vendor and product selection that have an impact on quality.

### SECTION 10. COMPLAINTS

#### 10.1 Overview

The laboratory considers an effective client complaint handling processes to be of significant business and strategic value. Listening to and documenting client concerns captures 'client knowledge' that enables our operations to continually improve processes and client satisfaction. An effective client complaint handling process also provides assurance to the data user that the laboratory will stand behind its data, service obligations and products.

A client complaint is any expression of dissatisfaction with any aspect of our business services (e.g., communications, responsiveness, data, reports, invoicing and other functions) expressed by any party, whether received verbally or in written form. Client inquiries, complaints or noted discrepancies are documented, communicated to management, and addressed promptly and thoroughly.

The laboratory has procedures for addressing both external and internal complaints with the goal of providing satisfactory resolution to complaints in a timely and professional manner.

The nature of the complaint is identified, documented and investigated, and an appropriate action is determined and taken. In cases where a client complaint indicates that an established policy or procedure was not followed, the QA Department must evaluate whether a special audit must be conducted to assist in resolving the issue. A written confirmation or letter to the client, outlining the issue and response taken is recommended as part of the overall action taken.

The process of complaint resolution and documentation utilizes the procedures outlined in Section 12 (Corrective Actions) and laboratory SOP TA-QA-0529 (Client Complaint Resolution).

#### 10.2 External Complaints

An employee that receives a complaint initiates the complaint resolution process by first documenting the complaint according to SOP TA-QA-0529.

Complaints fall into two categories: correctable and non-correctable. An example of a correctable complaint would be one where a report re-issue would resolve the complaint. An example of a non-correctable complaint would be one where a client complains that their data was repeatedly late. Non-correctable complaints should be reviewed for preventive action measures to reduce the likelihood of future occurrence and mitigation of client impact.

The general steps in the complaint handling process are:

- Receiving and Documenting Complaints

## Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 7  
Effective Date: 4/23/2014  
Page 47 of 165

- Complaint Investigation and Service Recovery
- Process Improvement

The laboratory shall inform the initiator of the complaint of the results of the investigation and the corrective action taken, if any.

### 10.3 Internal Complaints

Internal complaints include, but are not limited to: errors and non-conformances, training issues, internal audit findings, and deviations from methods. Corrective actions may be initiated by any staff member who observes a nonconformance and shall follow the procedures outlined in Section 12. In addition, Corporate Management, Sales and Marketing and IT may initiate a complaint by contacting the laboratory or through the corrective action system described in Section 12.

### 10.4 Management Review

The number and nature of client complaints is reported by the QA Manager to the laboratory and QA Director in the QA Monthly report. Monitoring and addressing the overall level and nature of client complaints and the effectiveness of the solutions is part of the Annual Management Review (Section 16).

## SECTION 11. CONTROL OF NON-CONFORMING WORK

### 11.1 Overview

When data discrepancies are discovered or deviations and departures from laboratory SOPs, policies and/or client requests have occurred, corrective action is taken immediately. First, the laboratory evaluates the significance of the nonconforming work. Then, a corrective action plan is initiated based on the outcome of the evaluation. If it is determined that the nonconforming work is an isolated incident, the plan could be as simple as adding a qualifier to the final results and/or making a notation in the case narrative. If it is determined that the nonconforming work is a systematic or improper practices issue, the corrective action plan could include a more in depth investigation and a possible suspension of an analytical method. In all cases, the actions taken are documented using the laboratory's corrective action system (refer to Section 12).

Due to the frequently unique nature of environmental samples, sometimes departures from documented policies and procedures are needed. When an analyst encounters such a situation, the problem is presented to the supervisor for resolution. The supervisor may elect to discuss it with the Technical Manager or have a representative contact the client to decide on a logical course of action. Once an approach is agreed upon, the analyst documents it using the laboratory's corrective action system described in Section 12. This information can then be supplied to the client in the form of a footnote or a case narrative with the report.

Project Management may encounter situations where a client may request that a special procedure be applied to a sample that is not standard lab practice. Based on a technical evaluation, the lab may accept or opt to reject the request based on technical or ethical merit. An example might be the need to report a compound that the lab does not normally report. The

## Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 7  
Effective Date: 4/23/2014  
Page 48 of 165

lab would not have validated the method for this compound following the procedures in Section 19. The client may request that the compound be reported based only on the calibration. Such a request would need to be approved by the Technical Manager and QA Manager, documented and included in the project folder. Deviations **must** also be noted on the final report with a statement that the compound is not reported in compliance with TNI (or the analytical method) requirements and the reason. Data being reported to a non-NELAC state would need to note the change made to how the method is normally run.

### 11.2 Responsibilities and Authorities

TestAmerica's Corporate SOP entitled Internal Investigation of Potential Data Discrepancies and Determination for Data Recall (SOP No. CW-L-S-002), outlines the general procedures for the reporting and investigation of data discrepancies and alleged incidents of misconduct or violations of TestAmerica's data integrity policies as well as the policies and procedures related to the determination of the potential need to recall data.

Under certain circumstances, the Laboratory Director, a Technical Manager, or a member of the QA team may authorize departures from documented procedures or policies. The departures may be a result of procedural changes due to the nature of the sample; a one-time procedure for a client; QC failures with insufficient sample to reanalyze, etc.. In most cases, the client will be informed of the departure prior to the reporting of the data. Any departures must be well documented using the laboratory's corrective action procedures. This information may also be documented in logbooks and/or data review checklists as appropriate. Any impacted data must be referenced in a case narrative and/or flagged with an appropriate data qualifier.

Any misrepresentation or possible misrepresentation of analytical data discovered by any laboratory staff member must be reported to facility Senior Management within 24-hours. The Senior Management staff is comprised of the Laboratory Director, the QA Manager, and the Technical Managers. The reporting of issues involving alleged violations of the company's Data Integrity or Manual Integration procedures **must** be conveyed to an Ethics and Compliance Officer (ECO), Corporate Quality Director and the laboratory's Quality Director within 24 hours of discovery.

Whether an inaccurate result was reported due to calculation or quantitation errors, data entry errors, improper practices, or failure to follow SOPs, the data must be evaluated to determine the possible effect.

The Laboratory Director, QA Manager, ECOs, Corporate Quality, General Managers and the Quality Directors have the authority and responsibility to halt work, withhold final reports, or suspend an analysis for due cause as well as authorize the resumption of work.

### 11.3 Evaluation of Significance and Actions Taken

For each nonconforming issue reported, an evaluation of its significance and the level of management involvement needed is made. This includes reviewing its impact on the final data, whether or not it is an isolated or systematic issue, and how it relates to any special client requirements.

TestAmerica's Corporate Data Investigation & Recall Procedure (SOP No. CW-L-S-002) distinguishes between situations when it would be appropriate for laboratory management to make the decision on the need for client notification (written or verbal) and data recall (report revision) and when the decision must be made with the assistance of the ECO's and Corporate Management. Laboratory level decisions are documented and approved using the laboratory's standard nonconformance/corrective action reporting in lieu of the data recall determination form contained in TestAmerica's Corporate SOP No. CW-L-S-002.

#### 11.4 Prevention of NonConforming Work

If it is determined that the nonconforming work could recur, further corrective actions must be made following the laboratory's corrective action system. On a monthly basis, the QA Department evaluates non-conformances to determine if any nonconforming work has been repeated multiple times. If so, the laboratory's corrective action process may be followed.

#### 11.5 Method Suspension / Restriction (Stop Work Procedures)

In some cases, it may be necessary to suspend/restrict the use of a method or target compound which constitutes significant risk and/or liability to the laboratory. Suspension/restriction procedures can be initiated by any of the persons noted in Section 11.2, Paragraph 5.

Prior to suspension/restriction, confidentiality will be respected, and the problem with the required corrective and preventive action will be stated in writing and presented to the Laboratory Director and Operations Manager.

The Laboratory Director and/or Operations Manager shall arrange for the appropriate personnel to meet with the QA Manager as needed. This meeting shall be held to confirm that there is a problem, that suspension/restriction of the method is required and will be concluded with a discussion of the steps necessary to bring the method/target or test fully back on line. In some cases, that may not be necessary if all appropriate personnel have already agreed there is a problem and there is agreement on the steps needed to bring the method, target or test fully back on line.

The QA Manager will also initiate a corrective action report as described in Section 12 if one has not already been started. A copy of any meeting notes and agreed upon steps should be faxed or e-mailed by the laboratory to the appropriate General Manager and member of Corporate QA. This fax/e-mail acts as notification of the incident.

After suspension/restriction, the lab will hold all reports to clients pending review. No faxing, mailing or distributing through electronic means may occur. The report must not be posted for viewing on the internet. It is the responsibility of the Laboratory Director to hold all reporting and to notify all relevant laboratory personnel regarding the suspension/restriction (e.g., Project Management, Log-in, etc...). Clients will NOT generally be notified at this time. Analysis may proceed in some instances depending on the non-conformance issue.

Within 72 hours, the QA Manager will determine if compliance is now met and reports can be released, OR determine the plan of action to bring work into compliance, and release work. A team, with all principals involved (Laboratory Director, Operations Manager, Technical Manager,

QA Manager) can devise a start-up plan to cover all steps from client notification through compliance and release of reports. Project Management and the Directors of Client Services and Sales and Marketing must be notified if clients must be notified or if the suspension/restriction affects the laboratory's ability to accept work. The QA Manager must approve start-up or elimination of any restrictions after all corrective action is complete. This approval is given by final signature on the completed corrective action report.

## SECTION 12. CORRECTIVE ACTION

### 12.1 Overview

A major component of TestAmerica's Quality Assurance (QA) Program is the problem investigation and feedback mechanism designed to keep the laboratory staff informed on quality related issues and to provide insight to problem resolution. When nonconforming work or departures from policies and procedures in the quality system or technical operations are identified, the corrective action procedure provides a systematic approach to assess the issues, restore the laboratory's system integrity, and prevent reoccurrence. Corrective actions are documented using Non-Conformance Memoranda (NCM) and Corrective Action Reports (CAR) (refer to SOP No. TA-QA-601 and Figure 12-1).

### 12.2 General

Problems within the quality system or within analytical operations may be discovered in a variety of ways, such as QC sample failures, internal or external audits, proficiency testing (PT) performance, client complaints, staff observation, etc.

The purpose of a corrective action system is to:

- Identify non-conformance events and assign responsibility(s) for investigating.
- Resolve non-conformance events and assign responsibility for any required corrective action.
- Identify systematic problems before they become serious.
- Identify and track client complaints and provide resolution.

**12.2.1 Non-Conformance Memorandum (NCM)** - is used to document the following types of corrective actions:

- Deviations from an established procedure or SOP
- QC outside of limits (non-matrix related)
- Isolated reporting / calculation errors
- Client complaints
- Discrepancies in materials / goods received vs. manufacturer packing slips.

**12.2.2 Corrective Action Report (CAR)** - is used to document the following types of corrective actions:

- Questionable trends that are found in the review of NCMs.
- Issues found while reviewing NCMs that warrant further investigation.
- Internal and external audit findings.
- Failed or unacceptable PT results.
- Corrective actions that cross multiple departments in the laboratory.
- Systematic reporting / calculation errors
- Client complaints
- Data recall investigations
- Identified poor process or method performance trends
- Excessive revised reports
- Health and Safety violations

This will provide background documentation to enable root cause analysis and preventive action.

### 12.3 Closed Loop Corrective Action Process

Any employee in the company can initiate a corrective action. There are four main components to a closed-loop corrective action process once an issue has been identified: Cause Analysis, Selection and Implementation of Corrective Actions (both short and long term), Monitoring of the Corrective Actions, and Follow-up. (refer to SOP No. TA-QA-0610 for more detail)

#### 12.3.1 Cause Analysis

- Upon discovery of a non-conformance event, the event must be defined and documented. An NCM or CAR must be initiated, someone is assigned to investigate the issue and the event is investigated for cause. Table 12-1 provides some general guidelines on determining responsibility for assessment.
- The cause analysis step is the key to the process as a long term corrective action cannot be determined until the cause is determined.
- If the cause is not readily obvious, the Technical Manager, Laboratory Director, or QA Manager (or QA designee) is consulted.

#### 12.3.2 Selection and Implementation of Corrective Actions

- Where corrective action is needed, the laboratory shall identify potential corrective actions. The action(s) most likely to eliminate the problem and prevent recurrence are selected and implemented. Responsibility for implementation is assigned.
- Corrective actions shall be to a degree appropriate to the magnitude of the problem identified through the cause analysis.
- Whatever corrective action is determined to be appropriate, the laboratory shall document and implement the changes. The NCM or CAR is used for this documentation.

#### 12.3.3 Root Cause Analysis

Root Cause Analysis is a class of problem solving (investigative) methods aimed at identifying the basic or causal factor(s) that underlie variation in performance or the occurrence of a

significant failure. The root cause may be buried under seemingly innocuous events, many steps preceding the perceived failure. At first glance, the immediate response is typically directed at a symptom and not the cause. Typically, root cause analysis would be best with three or more incidents to triangulate a weakness.

Systematically analyze and document the Root Causes of the more significant problems that are reported. Identify, track, and implement the corrective actions required to reduce the likelihood of recurrence of significant incidents. Trend the Root Cause data from these incidents to identify Root Causes that, when corrected, can lead to dramatic improvements in performance by eliminating entire classes of problems.

Identify the one event associated with problem and ask why this event occurred. Brainstorm the root causes of failures; for example, by asking why events occurred or conditions existed; and then why the cause occurred 5 consecutive times until you arrive at the root cause. Also, for each of these sub events or causes, ask why it occurred. Repeat the process for the other events associated with the incident.

Root cause analysis does not mean the investigation is over. Look at technique, or other systems outside the normal indicators. Often, creative thinking will find root causes that ordinarily would be missed, thus curtailing issues that could plague the laboratory or operation. (refer to SOP No. CA-Q-S-009 for more detail)

#### 12.3.4 Monitoring of the Corrective Actions

- The Technical Manager and QA Manager are responsible to ensure that the corrective action taken was effective.
- Ineffective actions are documented and re-evaluated until acceptable resolution is achieved. Technical Managers are accountable to the Operations Manager and Laboratory Director to ensure final acceptable resolution is achieved and documented appropriately.
- Each NCM and CAR is entered into a database for tracking purposes and a monthly summary of all corrective actions is printed out for review to aid in ensuring that the corrective actions have taken effect.
- The QA Manager reviews monthly NCMs and CARs for trends. Highlights are included in the QA monthly report (refer to Section 16). If a significant trend develops that adversely affects quality, an audit of the area is performed and corrective action implemented.
- Any out-of-control situations that are not addressed acceptably at the laboratory level may be reported to the Corporate Quality Director by the QA Manager, indicating the nature of the out-of-control situation and problems encountered in solving the situation.

#### 12.3.5 Follow-up Audits

- Follow-up audits may be initiated by the QA Manager and shall be performed as soon as possible when the identification of a nonconformance casts doubt on the laboratory's compliance with its own policies and procedures, or on its compliance with state or federal requirements.
- These audits often follow the implementation of the corrective actions to verify effectiveness. An additional audit would only be necessary when a critical issue or risk to business is discovered.

(Also refer to Section 15.1.4, Special Audits.)

## 12.4 Technical Corrective Actions

In addition to providing acceptance criteria and specific protocols for technical corrective actions in the method SOPs, the laboratory has general procedures to be followed to determine when departures from the documented policies and procedures and quality control have occurred (refer to Section 11). The documentation of these procedures is through the use of an NCM.

Table 12-1 includes examples of general technical corrective actions. For specific criteria and corrective actions, refer to the analytical methods or specific method SOPs. The laboratory may also maintain Work Instructions on these items that are available upon request.

Table 12-1 provides some general guidelines for identifying the individual(s) responsible for assessing each QC type and initiating corrective action. The table also provides general guidance on how a data set should be treated if associated QC measurements are unacceptable. Specific procedures are included in Method SOPs, Work Instructions and QAM Sections 19 and 20. All corrective actions are reviewed monthly, at a minimum, by the QA Manager and highlights are included in the QA monthly report.

To the extent possible, samples shall be reported only if all quality control measures are acceptable. If the deficiency does not impair the usability of the results, data will be reported with an appropriate data qualifier and/or the deficiency will be noted in the case narrative. Where sample results may be impaired, the Project Manager is notified by an NCM and appropriate corrective action (e.g., reanalysis) is taken and documented.

## 12.5 Basic Corrections

When mistakes occur in records, each mistake shall be crossed-out, [not obliterated (e.g. no white-out)], and the correct value entered alongside. All such corrections shall be initialed (or signed) and dated by the person making the correction. In the case of records stored electronically, the original "uncorrected" file must be maintained intact and a second "corrected" file is created.

This same process applies to adding additional information to a record. All additions made later than the initial must also be initialed (or signed) and dated.

When corrections are due to reasons other than obvious transcription errors, the reason for the corrections (or additions) shall also be documented.

Figure 12-1.  
Example - Non Conformance Memo

01/09/2008	<b>NonConformance Memo Report</b>	Page 1 of 1
NCM ID: 580-5190		Date Opened: 11/30/2007
NCM Type: CCV - %D, High, Sample ND		Date Closed:
Lab Section: Pesticides and PCBs		Date Verified:
<b>Narrative</b> The closing continuing calibration verification (CCV) for Endosulfan II recovered above the upper control limit. The samples associated with this CCV were non-detects for the affected analytes; therefore, the data have been reported.		
<b>Affected Items</b>		
Description	Project Manager	
580-8247-E-1-A	Curbow, Heather	
LCS 580-26165/2-A	Curbow, Heather	
580-8247-C-2-A	Curbow, Heather	
580-8247-E-1-A	Curbow, Heather	
LCS 580-26165/2-A	Curbow, Heather	
MB 580-26165/1-A	Curbow, Heather	
LCSD 580-26165/3-A	Curbow, Heather	
580-8247-C-2-A	Curbow, Heather	
MB 580-26165/1-A	Curbow, Heather	
LCSD 580-26165/3-A	Curbow, Heather	
<b>Notifications</b>		
User Full Name	Date Received	Date
Mott, Christina M	11/30/2007	
	Notice Level	Verification Type
	Level 1	Review
<b>Job Case Narratives</b>		
Job: 580-J8247-1		
Job Narrative		
580-J8247-1		
Comments		
No additional comments.		
Receipt		
All samples were received in good condition within temperature requirements.		
GC/MS VOA		
No analytical or quality issues were noted.		
GC/MS Semi VOA		
No analytical or quality issues were noted.		
GC Semi VOA		
Method 8081: The RPD between the results of the two columns for heptachlor and DDT was >40%. This is due to matrix effect; the lower of the two values (presumably free from matrix interference) is reported as primary.		
No other analytical or quality issues were noted.		
Metals		
No analytical or quality issues were noted.		
General Chemistry		
No analytical or quality issues were noted.		
Organic Prep		
No analytical or quality issues were noted.		



Table 12-1. Example – General Corrective Action Procedures

QC Activity (Individual Responsible for Initiation/Assessment)	Acceptance Criteria	Recommended Corrective Action
Initial Instrument Blank  (Analyst)	- Instrument response < MDL.	- Prepare another blank. - If same response, determine cause of contamination: reagents, environment, instrument equipment failure, etc..
Initial Calibration Standards  (Analyst, Technical Manager(s))	- Correlation coefficient > 0.99 or standard concentration value. - % Recovery within acceptance range. - See details in Method SOP.	- Reanalyze standards. - If still unacceptable, remake standards and recalibrate instrument.
Independent Calibration Verification (Second Source)  (Analyst, Technical Manager(s))	- % Recovery within control limits.	- Remake and reanalyze standard. - If still unacceptable, then remake calibration standards or use new primary standards and recalibrate instrument.
Continuing Calibration Standards  (Analyst, Data Reviewer)	% Recovery within control limits.	- Reanalyze standard. - If still unacceptable, then recalibrate and rerun affected samples.
Matrix Spike / Matrix Spike Duplicate (MS/MSD)  (Analyst, Data Reviewer)	- % Recovery within limits documented in LIMS.	- If the acceptance criteria for duplicates or matrix spikes are not met because of matrix interferences, the acceptance of the analytical batch is determined by the validity of the LCS. - If the LCS is within acceptable limits the batch is acceptable. - The results of the duplicates, matrix spikes and the LCS are reported with the data set. - For matrix spike or duplicate results outside criteria, the data for that sample shall be reported with qualifiers.

QC Activity (Individual Responsible for Initiation/Assessment)	Acceptance Criteria	Recommended Corrective Action
Laboratory Control Sample (LCS)  (Analyst, Data Reviewer)	- % Recovery within limits specified in LIMS.	- Batch must be re-prepared and re-analyzed. This includes any allowable marginal exceedance. When not using marginal exceedances, the following exceptions apply: 1) when the acceptance criteria for the positive control are exceeded high (i.e., high bias) and there are associated samples that are non-detects, then those non-detects may be reported with data qualifying codes; 2) when the acceptance criteria for the positive control are exceeded low (i.e., low bias), those sample results may be reported if they exceed a maximum regulatory limit/decision level with data qualifying codes.  <b>Note:</b> If there is insufficient sample or the holding time cannot be met, contact client and report with flags.
Surrogates  (Analyst, Data Reviewer)	- % Recovery within limits of method or within three standard deviations of the historical mean.	- Individual sample must be repeated. Place comment in LIMS. - Surrogate results outside criteria shall be reported with qualifiers.
Method Blank (MB)  (Analyst, Data Reviewer)	< Reporting Limit <sup>1</sup>	- Reanalyze blank. - If still positive, determine source of contamination. If necessary, reprocess (i.e. digest or extract) entire sample batch. Report blank results. - Qualify the result(s) if the concentration of a targeted analyte in the MB is at or above the reporting limit AND is > 1/10 of the amount measured in the sample.
Proficiency Testing (PT) Samples  (QA Manager, Technical Manager(s))	- Criteria supplied by PT Supplier.	- Any failures or warnings must be investigated for cause. Failures may result in the need to repeat a PT sample to show the problem is corrected.
Internal / External Audits  (QA Manager, Technical Manager(s), Laboratory Director)	- Defined in Quality System documentation such as SOPs, QAM, etc..	- Non-conformances must be investigated through CAR system and necessary corrections must be made.



QC Activity (Individual Responsible for Initiation/Assessment)	Acceptance Criteria	Recommended Corrective Action
Reporting / Calculation Errors  (Depends on issue – possible individuals include: Analysts, Data Reviewers, Project Managers, Technical Managers, QA Manager, Corporate QA, Corporate Management)	- SOP CW-L-S-002, Internal Investigation of Potential Data Discrepancies and Determination for Data Recall.	- Corrective action is determined by type of error. Follow the procedures in SOP CW-L-S-002 or Lab SOP TA-QA- 0610 Laboratory Corrective Action Procedures.
Client Complaints  (Project Managers, Lab Director, Sales and Marketing)	-	- Corrective action is determined by the type of complaint. For example, a complaint regarding an incorrect address on a report will result in the report being corrected and then follow- up must be performed on the reasons the address was incorrect (e.g., database needs to be updated).
QA Monthly Report (Refer to Section 16 for an example)  (QA Manager, Lab Director, Technical Manager(s))	- QAM, SOPs.	- Corrective action is determined by the type of issue. For example, CARs for the month are reviewed and possible trends are investigated.
Health and Safety Violation  (Safety Officer, Lab Director, Technical Manager(s))	- Environmental Health and Safety (EHS) Manual.	- Non-conformance is investigated and corrected through CAR system.

**Note:**

1. Except as noted below for certain compounds, the method blank should be below the detection limit or for DoD projects the method blank should be below  $\frac{1}{2}$  the RL. Concentrations up to five times the reporting limit will be allowed for the ubiquitous laboratory and reagent contaminants: methylene chloride, toluene, acetone, 2-butanone and phthalates provided they appear in similar levels in the reagent blank and samples. This allowance presumes that the detection limit is significantly below any regulatory limit to which the data are to be compared and that blank subtraction will not occur. For benzene and ethylene dibromide (EDB) and other analytes for which regulatory limits are extremely close to the detection limit, the method blank must be below the method detection limit.

**SECTION 13. PREVENTIVE ACTION / IMPROVEMENT****13.1 Overview**

The laboratory's preventive action programs improve or eliminate potential causes of nonconforming product and/or nonconformance to the quality system. This preventive action process is a proactive and continuous process of improvement activities that can be initiated through feedback from clients, employees, business providers, and affiliates. The QA Department has the overall responsibility to ensure that the preventive action process is in place, and that relevant information on actions is submitted for management review.

Dedicating resources to an effective preventive action system emphasizes the laboratory's commitment to its Quality Program. It is beneficial to identify and address negative trends before they develop into complaints, problems and corrective actions. Additionally, customer service and client satisfaction can be improved through continuous improvements to laboratory systems.

Opportunities for improvement may be discovered during management reviews, the monthly QA Metrics Report, evaluation of internal or external audits, results & evaluation of proficiency testing (PT) performance, data analysis & review processing operations, client complaints, staff observation, etc..

The monthly QA Metrics Report shows performance indicators in all areas of the laboratory and quality system. These areas include revised reports, corrective actions, audit findings, internal auditing and data authenticity audits, client complaints, PT samples, holding time violations, SOPs, ethics training, etc.. These metrics are used in evaluating the management and quality system performance on an ongoing basis and provide a tool for identifying areas for improvement.

The laboratory's corrective action process is integral to implementation of preventive actions. A critical piece of the corrective action process is the implementation of actions to prevent further occurrence of a non-compliance event. Historical review of corrective action provides a valuable mechanism for identifying preventive action opportunities.

**13.1.1 The following elements are part of a preventive action system:**

- Identification of an opportunity for preventive action.
- Process for the preventive action.
- Define the measurements of the effectiveness of the process once undertaken.
- Execution of the preventive action.
- Evaluation of the plan using the defined measurements.
- Verification of the effectiveness of the preventive action.
- Close-Out by documenting any permanent changes to the Quality System as a result of the Preventive Action. Documentation of Preventive Action is incorporated into the monthly QA reports, corrective action process and management review.

13.1.2 Any Preventive Actions undertaken or attempted shall be taken into account during the annual Management Systems Review (Section 16). A highly detailed report is not required; however, a summary of successes and failures within the preventive action program is sufficient to provide management with a measurement for evaluation.

### 13.2 Management of Change

The Management of Change process is designed to manage significant events and changes that occur within the laboratory. Through these procedures, the potential risks inherent with a new event or change are identified and evaluated. The risks are minimized or eliminated through pre-planning and the development of preventive measures. The types of changes covered under this system include: Facility Changes, Major Accreditation Changes, Addition or Deletion to Division's Capabilities or Instrumentation, Key Personnel Changes, Laboratory Information Management System (LIMS) changes. This process is discussed in further detail in SOP TA-QA-0530, Management of Change.

## SECTION 14. CONTROL OF RECORDS

The laboratory maintains a records management system appropriate to its needs and that complies with applicable standards or regulations as required. The system produces unequivocal, accurate records that document all laboratory activities. The laboratory retains all original observations, calculations and derived data, calibration records and a copy of the analytical report for a minimum of five years after it has been issued.

### 14.1 Overview

The laboratory has established procedures for identification, collection, indexing, access, filing, storage, maintenance and disposal of quality and technical records. A record index is listed in Table 14-1. Quality records are maintained by the QA department in a database, which is backed up as part of the regular laboratory backup. Records are of two types; either electronic or hard copy paper formats depending on whether the record is computer or hand generated (some records may be in both formats, in this instance the electronic copy will be considered to be the vital record). Technical records are maintained by the Technical Managers.

Table 14-1. Record Index<sup>1</sup>

	Record Types <sup>1</sup> :	Retention Time:
<b>Technical Records</b>	<ul style="list-style-type: none"> <li>- Raw Data</li> <li>- Logbooks<sup>2</sup></li> <li>- Standards</li> <li>- Certificates</li> <li>- Analytical Records</li> <li>- MDLs/IDLs/DOCs</li> <li>- Lab Reports</li> </ul>	5 Years from analytical report issue*

	Record Types <sup>1</sup> :	Retention Time:
<b>Official Documents</b>	<ul style="list-style-type: none"> <li>- Quality Assurance Manual (QAM)</li> <li>- Work Instructions</li> <li>- Policies</li> <li>- SOPs</li> <li>- Policy Memorandums</li> <li>- Manuals</li> </ul>	5 Years from document retirement date*
<b>QA Records</b>	<ul style="list-style-type: none"> <li>- Internal &amp; External Audits/Responses</li> <li>- Certifications</li> <li>- Corrective/Preventive Actions</li> <li>- Management Reviews</li> <li>- Method &amp; Software Validation / Verification Data</li> <li>- Data Investigation</li> </ul>	5 Years from archival*  <b>Data Investigation:</b> 5 years or the life of the affected raw data storage whichever is greater (beyond 5 years if ongoing project or pending investigation)
<b>Project Records</b>	<ul style="list-style-type: none"> <li>- Sample Receipt &amp; COC Documentation</li> <li>- Contracts and Amendments</li> <li>- Correspondence</li> <li>- QAPP</li> <li>- SAP</li> <li>- Telephone Logbooks</li> <li>- Lab Reports</li> </ul>	5 Years from analytical report issue*
<b>Administrative Records</b>	Finance and Accounting	10 years
	EH&S Manual, Permits	7 years
	<b>Disposal Records, Permits</b>	<b>Indefinitely</b>
	Employee Handbook	Indefinitely
	Personnel files, Employee Signature & Initials, Administrative Training Records (e.g., Ethics)	Refer to HR Manual
	Administrative Policies Technical Training Records	7 years

<sup>1</sup> Record Types encompass hardcopy and electronic records.

<sup>2</sup> Examples of Logbook types: Maintenance, Instrument Run, Preparation (standard and samples), Standard and Reagent Receipt, Archiving, Balance Calibration, Temperature (hardcopy or electronic records).

\* Exceptions listed in Table 14-2.

14.1.1 All records are stored and retained in such a way that they are secure and readily retrievable at the laboratory or the Iron Mountain data storage facility that provide a suitable environment to prevent damage or deterioration and to prevent loss. All records shall be protected against fire, theft, loss, environmental deterioration, and vermin. In the case of electronic records, electronic or magnetic sources, storage media are protected from deterioration caused by magnetic fields and/or electronic deterioration.

Access to the data is limited to laboratory and company employees and shall be documented with an access log. Records archived off-site are stored in a secure location where a record is maintained of any entry into the storage facility. Whether on-site or off-site storage is used, logs are maintained in each storage box to note removal and return of records. Retention of records are maintained on-site at the laboratory for at least 6 months after their generation and moved

offsite for the remainder of the required storage time. Records are maintained for a minimum of five years unless otherwise specified by a client or regulatory requirement.

For raw data and project records, record retention shall be calculated from the date the project report is issued. For other records, such as Controlled Documents, QA, or Administrative Records, the retention time is calculated from the date the record is formally retired. Records related to the programs listed in Table 14-2 have lengthier retention requirements and are subject to the requirements in Section 14.1.3.

#### 14.1.2 Programs with Longer Retention Requirements

Some regulatory programs have longer record retention requirements than the standard record retention time. These are detailed in Table 14-2 with their retention requirements. In these cases, the longer retention requirement is enacted. If special instructions exist such that client data cannot be destroyed prior to notification of the client, the container or box containing that data is marked as to who to contact for authorization prior to destroying the data.

**Table 14-2. Example: Special Record Retention Requirements**

Program	Retention Requirement
Drinking Water – All States	5 years (project records) 10 years - Radiochemistry (project records)
Drinking Water Lead and Copper Rule	12 years (project records)
Alaska	10 years
Louisiana – All	10 years
Navy Facilities Engineering Service Center (NFESC)	10 years
TSCA - 40 CFR Part 792	10 years after publication of final test rule or negotiated test agreement

<sup>1</sup>Note: Extended retention requirements must be noted with the archive documents or addressed in facility-specific records retention procedures.

**14.1.3** The laboratory has procedures to protect and back-up records stored electronically and to prevent unauthorized access to or amendment of these records. All analytical data is maintained as hard copy or in a secure readable electronic format. For analytical reports that are maintained as copies in PDF format, refer to Section 19.14.1 for more information. Additional information can also be reference in SOP TA-QA-0506 Archiving Data and Reports.

**14.1.4** The record keeping system allows for historical reconstruction of all laboratory activities that produced the analytical data, as well as rapid recovery of historical data. The history of the sample from when the laboratory took possession of the samples must be readily understood through the documentation. This shall include inter-laboratory transfers of samples and/or extracts.

- The records include the identity of personnel involved in sampling, sample receipt, preparation, or testing. All analytical work contains the initials (at least) of the personnel involved. The laboratory's copy of the COC is stored with the invoice and the work order

sheet generated by the LIMS in a job folder. The chain of custody would indicate the name of the sampler. If any sampling notes are provided with a work order, they are kept with this package.

- All information relating to the laboratory facilities equipment, analytical test methods, and related laboratory activities, such as sample receipt, sample preparation, or data verification are documented.
- The record keeping system facilitates the retrieval of all working files and archived records for inspection and verification purposes (e.g., set format for naming electronic files, set format for what is included with a given analytical data set as per SOP TA-QA-0506 Archiving Reports and Report File Maintenance). Instrument data is stored sequentially by instrument. A given day's analyses are maintained in the order of the analysis. Run logs are maintained for each instrument or method; a copy of each day's run log or instrument sequence is stored with the data to aid in re-constructing an analytical sequence. Where an analysis is performed without an instrument, bound logbooks or bench sheets are used to record and file data. Standard and reagent information is recorded in logbooks or entered into the LIMS for each method as required.
- Changes to hardcopy records shall follow the procedures outlined in Section 12 and 19. Changes to electronic records in LIMS or instrument data are recorded in audit trails.
- The reason for a signature or initials on a document is clearly indicated in the records such as "sampled by," "prepared by," "reviewed by", or "analyzed by".
- All generated data except those that are generated by automated data collection systems, are recorded directly, promptly and legibly in permanent dark ink.
- Hard copy data may be scanned into PDF format for record storage as long as the scanning process can be verified in order to ensure that no data is lost and the data files and storage media must be tested to verify the laboratory's ability to retrieve the information prior to the destruction of the hard copy that was scanned. The procedure for this verification can be found in SOP TA-QA-0506 Archiving Reports and Report File Maintenance.
- Also refer to Section 19.14.1 'Computer and Electronic Data Related Requirements'.

#### 14.2 Technical and Analytical Records

**14.2.1** The laboratory retains records of original observations, derived data and sufficient information to establish an audit trail, calibration records, staff records and a copy of each analytical report issued, for a minimum of five years unless otherwise specified by a client or regulatory requirement. The records for each analysis shall contain sufficient information to enable the analysis to be repeated under conditions as close as possible to the original. The records shall include the identity of laboratory personnel responsible for the performance of each analysis and reviewing results.

**14.2.2** Observations, data and calculations are recorded real-time and are identifiable to the specific task.

## Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 7  
Effective Date: 4/23/2014  
Page 63 of 165

**14.2.3** Changes to hardcopy records shall follow the procedures outlined in Section 12 and 19. Changes to electronic records in LIMS or instrument data are recorded in audit trails.

The essential information to be associated with analysis, such as strip charts, tabular printouts, computer data files, analytical notebooks, and run logs, include:

- laboratory sample ID code;
- Date of analysis; Time of Analysis is also required if the holding time is seventy-two (72) hours or less, or when time critical steps are included in the analysis (e.g., drying times, incubations, etc.); instrumental analyses have the date and time of analysis recorded as part of their general operations. Where a time critical step exists in an analysis, location for such a time is included as part of the documentation in a specific logbook or on a benchsheet or in the LIMS.
- Instrumentation identification and instrument operating conditions/parameters. Operating conditions/parameters are typically recorded in instrument maintenance logs where available.
- analysis type;
- all manual calculations and manual integrations;
- analyst's or operator's initials/signature;
- sample preparation including cleanup, separation protocols, incubation periods or subculture, ID codes, volumes, weights, instrument printouts, meter readings, calculations, reagents;
- test results;
- standard and reagent origin, receipt, preparation, and use;
- calibration criteria, frequency and acceptance criteria;
- data and statistical calculations, review, confirmation, interpretation, assessment and reporting conventions;
- quality control protocols and assessment;
- electronic data security, software documentation and verification, software and hardware audits, backups, and records of any changes to automated data entries; and
- Method performance criteria including expected quality control requirements. These are indicated both in the LIMS and on specific analytical report formats.

### **14.3 Laboratory Support Activities**

In addition to documenting all the above-mentioned activities, the following are retained QA records and project records (previous discussions in this section relate where and how these data are stored):

- all original raw data, whether hard copy or electronic, for calibrations, samples and quality control measures, including analysts' work sheets and data output records (chromatograms, strip charts, and other instrument response readout records);
- a written description or reference to the specific test method used which includes a

## Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 7  
Effective Date: 4/23/2014  
Page 64 of 165

description of the specific computational steps used to translate parametric observations into a reportable analytical value;

- copies of final reports;
- archived SOPs;
- correspondence relating to laboratory activities for a specific project;
- all corrective action reports, audits and audit responses;
- proficiency test results and raw data; and
- results of data review, verification, and crosschecking procedures.

### **14.3.1 Sample Handling Records**

Records of all procedures to which a sample is subjected while in the possession of the laboratory are maintained. These include but are not limited to records pertaining to:

- sample preservation including appropriateness of sample container and compliance with holding time requirement;
- sample identification, receipt, acceptance or rejection and login;
- sample storage and tracking including shipping receipts, sample transmittal / COC forms; and
- procedures for the receipt and retention of samples, including all provisions necessary to protect the integrity of samples.

### **14.4 Administrative Records**

The laboratory also maintains the administrative records in either electronic or hard copy form. Refer to Table 14-1.

### **14.5 Records Management, Storage and Disposal**

All records (including those pertaining to test equipment), certificates and reports are safely stored, held secure and in confidence to the client. Certification related records are available upon request.

All information necessary for the historical reconstruction of data is maintained by the laboratory. Records that are stored only on electronic media must be supported by the hardware and software necessary for their retrieval.

Records that are stored or generated by computers or personal computers have hard copy, write-protected backup copies, or an electronic audit trail controlling access.

The laboratory has a record management system (a.k.a., document control) for control of laboratory notebooks, instrument logbooks, standards logbooks, and records for data reduction, validation, storage and reporting. Laboratory notebooks are issued on a per analysis basis, and are numbered sequentially. All data are recorded sequentially within a series of sequential

notebooks. Bench sheets are filed sequentially. Standards are maintained in the LIMS – no logbooks are used to record that data. Records are considered archived when noted as such in the records management system (a.k.a., document control.)

#### 14.5.1 Transfer of Ownership

In the event that the laboratory transfers ownership or goes out of business, the laboratory shall ensure that the records are maintained or transferred according to client's instructions. Upon ownership transfer, record retention requirements shall be addressed in the ownership transfer agreement and the responsibility for maintaining archives is clearly established. In addition, in cases of bankruptcy, appropriate regulatory and state legal requirements concerning laboratory records must be followed. In the event of the closure of the laboratory, all records will revert to the control of the corporate headquarters. Should the entire company cease to exist, as much notice as possible will be given to clients and the accrediting bodies who have worked with the laboratory during the previous 5 years of such action.

#### 14.5.2 Records Disposal

Records are removed from the archive and destroyed after 5 years unless otherwise specified by a client or regulatory requirement. On a project specific or program basis, clients may need to be notified prior to record destruction. Records are destroyed in a manner that ensures their confidentiality such as shredding, mutilation or incineration. (Refer to Tables 14-1 and 14-2).

Electronic copies of records must be destroyed by erasure or physically damaging off-line storage media so no records can be read.

If a third party records management company is hired to dispose of records, a "Certificate of Destruction" is required.

### SECTION 15. AUDITS

#### 15.1 Internal Audits

Internal audits are performed to verify that laboratory operations comply with the requirements of the lab's quality system and with the external quality programs under which the laboratory operates. Audits are planned and organized by the QA staff. Personnel conducting the audits should be independent of the area being evaluated. Auditors will have sufficient authority, access to work areas, and organizational freedom necessary to observe all activities affecting quality and to report the assessments to laboratory management and, when requested, to corporate management.

Audits are conducted and documented as described in the TestAmerica Corporate SOP on performing Internal Auditing, SOP No. CA-Q-S-004. The types and frequency of routine internal audits are described in Table 15-1. Special or ad hoc assessments may be conducted as needed under the direction of the QA staff.

**Table 15-1. Types of Internal Audits and Frequency**

Description	Performed by	Frequency
Quality Systems Audits	QA Department, QA approved designee, or Corporate QA	All areas of the laboratory annually
Method Audits	Joint responsibility: a) QA Manager or designee b) Technical Manager or Designee (Refer to CA-Q-S-004)	Methods Audits Frequency: <b>100% of methods annually</b>
Special	QA Department or Designee	Surveillance or spot checks performed as needed, e.g., to confirm corrective actions from other audits.
Performance Testing	Analysts with QA oversight	Two successful per year for each NELAC field of testing or as dictated by regulatory requirements

#### 15.1.1 Annual Quality Systems Audit

An annual quality systems audit is required to ensure compliance to analytical methods and SOPs, TestAmerica's Data Integrity and Ethics Policies, NELAC quality systems, client and state requirements, and the effectiveness of the internal controls of the analytical process, including but not limited to data review, quality controls, preventive action and corrective action. The completeness of earlier corrective actions is assessed for effectiveness & sustainability. The audit is divided into sections for each operating or support area of the lab, and each section is comprehensive for a given area. The area audits may be performed on a rotating schedule throughout the year to ensure adequate coverage of all areas. This schedule may change as situations in the laboratory warrant.

#### 15.1.2 Method Audits

Method audits are based on client projects, associated sample delivery groups, and the methods performed. Reported results are compared to raw data to verify the authenticity of results. The validity of calibrations and QC results are compared to data qualifiers, footnotes, and case narratives. Documentation is assessed by examining run logs and records of manual integrations. Manual calculations are checked. Where possible, electronic audit miner programs (e.g., MintMiner and Chrom AuditMiner) are used to identify unusual manipulations of the data deserving closer scrutiny. Method audits will include all methods annually.

#### 15.1.3 SOP Method Compliance

Compliance of all SOPs with the source methods and compliance of the operational groups with the SOPs will be assessed by the Technical Manager or qualified designee annually. It is also recommended that the work of each newly hired analyst is assessed within 3 months of working independently, (e.g., completion of method IDOC). In addition, as analysts add methods to their capabilities, (new IDOC) reviews of the analyst work products will be performed within 3 months of completing the documented training.

**15.1.4 Special Audits**

Special audits are conducted on an as needed basis, generally as a follow up to specific issues such as client complaints, corrective actions, PT results, data audits, system audits, validation comments, regulatory audits or suspected ethical improprieties. Special audits are focused on a specific issue, and report format, distribution, and timeframes are designed to address the nature of the issue.

**15.1.5 Performance Testing**

The laboratory participates semi-annually in performance audits conducted through the analysis of PT samples provided by a third party. The laboratory generally participates in the following types of PT studies: Drinking Water, Non Potable Water, and Soil.

It is TestAmerica's policy that PT samples be treated as typical samples in the production process. Furthermore, where PT samples present special or unique problems, in the regular production process they may need to be treated differently, as would any special or unique request submitted by any client. The QA Manager must be consulted and in agreement with any decisions made to treat a PT sample differently due to some special circumstance.

Written responses to unacceptable PT results are required. In some cases it may be necessary for blind QC samples to be submitted to the laboratory to show a return to control.

**15.2 External Audits**

External audits are performed when certifying agencies or clients conduct on-site inspections or submit performance testing samples for analysis. It is TestAmerica's policy to cooperate fully with regulatory authorities and clients. The laboratory makes every effort to provide the auditors with access to personnel, documentation, and assistance. Laboratory supervisors are responsible for providing corrective actions to the QA Manager who coordinates the response for any deficiencies discovered during an external audit. Audit responses are due in the time allotted by the client or agency performing the audit. When requested, a copy of the audit report and the labs corrective action plan will be forwarded to Corporate Quality.

The laboratory cooperates with clients and their representatives to monitor the laboratory's performance in relation to work performed for the client. The client may only view data and systems related directly to the client's work. All efforts are made to keep other client information confidential.

**15.2.1 Confidential Business Information (CBI) Considerations**

During on-site audits, auditors may come into possession of information claimed as business confidential. A business confidentiality claim is defined as "a claim or allegation that business information is entitled to confidential treatment for reasons of business confidentiality or a request for a determination that such information is entitled to such treatment." When information is claimed as business confidential, the laboratory must place on (or attach to) the information at the time it is submitted to the auditor, a cover sheet, stamped or typed legend or other suitable form of notice, employing language such as "trade secret", "proprietary" or "company confidential". Confidential portions of documents otherwise non-confidential must be clearly identified. CBI may be purged of references to client identity by the responsible

laboratory official at the time of removal from the laboratory. However, sample identifiers may not be obscured from the information. Additional information regarding CBI can be found in within the 2009 TNI standards.

**15.3 Audit Findings**

Audit findings are documented using the corrective action process and database. The laboratory's corrective action responses for both types of audits may include action plans that could not be completed within a predefined timeframe. In these instances, a completion date must be set and agreed to by operations management and the QA Manager.

Developing and implementing corrective actions to findings is the responsibility of the Technical Manager where the finding originated. Findings that are not corrected by specified due dates are reported monthly to management in the QA monthly report. When requested, a copy of the audit report and the labs corrective action plan will be forwarded to Corporate Quality.

If any audit finding casts doubt on the effectiveness of the operations or on the correctness or validity of the laboratory's test results, the laboratory shall take timely corrective action, and shall notify clients in writing if the investigations show that the laboratory results have been affected. Once corrective action is implemented, a follow-up audit is scheduled to ensure that the problem has been corrected.

Clients must be notified promptly in writing, of any event such as the identification of defective measuring or test equipment that casts doubt on the validity of results given in any test report or amendment to a test report. The investigation must begin within 24-hours of discovery of the problem and all efforts are made to notify the client within two weeks after the completion of the investigation.

**SECTION 16. MANAGEMENT REVIEWS****16.1 Quality Assurance Report**

A comprehensive QA Report shall be prepared each month by the laboratory's QA Department and forwarded to the Laboratory Director, Technical Managers, Operation Manager, their Quality Director as well as the General Manager. All aspects of the QA system are reviewed to evaluate the suitability of policies and procedures. During the course of the year, the Laboratory Director, General Manager or Corporate QA may request that additional information be added to the report.

On a monthly basis, Corporate QA compiles information from all the monthly laboratory reports. The Corporate Quality Directors prepare a report that includes a compilation of all metrics and notable information and concerns regarding the QA programs within the laboratories. The report also includes a listing of new regulations that may potentially impact the laboratories. This report is presented to the Senior Management Team and General Managers.



## 16.2 Annual Management Review

The senior lab management team (Laboratory Director, Technical Managers, Operation Manager, QA Manager) conducts a review annually of its quality systems and LIMS to ensure its continuing suitability and effectiveness in meeting client and regulatory requirements and to introduce any necessary changes or improvements. It will also provide a platform for defining goals, objectives and action items that feed into the laboratory planning system. Corporate Operations and Corporate QA personnel is included in this meeting at the discretion of the Laboratory Director. The LIMS review consists of examining any audits, complaints or concerns that have been raised through the year that are related to the LIMS. The laboratory will summarize any critical findings that can not be solved by the lab and report them to Corporate IT.

This management systems review (Corporate SOP No. CA-Q-S-008 & Work Instruction No. CA-Q-WI-020) uses information generated during the preceding year to assess the "big picture" by ensuring that routine actions taken and reviewed on a monthly basis are not components of larger systematic concerns. The monthly review should keep the quality systems current and effective, therefore, the annual review is a formal senior management process to review specific existing documentation. Significant issues from the following documentation are compiled or summarized by the QA Manager prior to the review meeting:

- Matters arising from the previous annual review.
- Prior Monthly QA Reports issues.
- Laboratory QA Metrics.
- Review of report reissue requests.
- Review of client feedback and complaints.
- Issues arising from any prior management or staff meetings.
- Minutes from prior senior lab management meetings. Issues that may be raised from these meetings include:
  - Adequacy of staff, equipment and facility resources.
  - Adequacy of policies and procedures.
  - Future plans for resources and testing capability and capacity.
- The annual internal double blind PT program sample performance (if performed),
- Compliance to the Ethics Policy and Data Integrity Plan. Including any evidence/incidents of inappropriate actions or vulnerabilities related to data Integrity.

A report is generated by the QA Manager and management. The report is distributed to the appropriate General Manager and the Quality Director. The report includes, but is not limited to:

- The date of the review and the names and titles of participants.
- A reference to the existing data quality related documents and topics that were reviewed.
- Quality system or operational changes or improvements that will be made as a result of the review [e.g., an implementation schedule including assigned responsibilities for the changes (Action Table)].

Changes to the quality systems requiring update to the laboratory QA Manual shall be included in the next revision of the QA Manual.

## 16.3 Potential Integrity Related Managerial Reviews

Potential integrity issues (data or business related) must be handled and reviewed in a confidential manner until such time as a follow-up evaluation, full investigation, or other appropriate actions have been completed and issues clarified. TestAmerica's Corporate Data Investigation/Recall SOP shall be followed (SOP No. CW-L-S-002). All investigations that result in finding of inappropriate activity are documented and include any disciplinary actions involved, corrective actions taken, and all appropriate notifications of clients.

TestAmerica's CEO, VP of Quality, Technical & Operations Support, General Managers and Quality Directors receive a monthly report from the Corporate Quality Director summarizing any current data integrity or data recall investigations. The General Managers are also made aware of progress on these issues for their specific labs.

## SECTION 17. PERSONNEL

### 17.1 Overview

The laboratory's management believes that its highly qualified and professional staff is the single most important aspect in assuring a high level of data quality and service. The staff consists of professionals and support personnel as outlined in the organization chart in Figure 4-1.

All personnel must demonstrate competence in the areas where they have responsibility. Any staff that is undergoing training shall have appropriate supervision until they have demonstrated their ability to perform their job function on their own. Staff shall be qualified for their tasks based on appropriate education, training, experience and/or demonstrated skills as required.

The laboratory employs sufficient personnel with the necessary education, training, technical knowledge and experience for their assigned responsibilities.

All personnel are responsible for complying with all QA/QC requirements that pertain to the laboratory and their area of responsibility. Each staff member must have a combination of experience and education to adequately demonstrate a specific knowledge of their particular area of responsibility. Technical staff must also have a general knowledge of lab operations, test methods, QA/QC procedures and records management.

Laboratory management is responsible for formulating goals for lab staff with respect to education, training and skills and ensuring that the laboratory has a policy and procedures for identifying training needs and providing training of personnel. The training shall be relevant to the present and anticipated responsibilities of the lab staff.

The laboratory only uses personnel that are employed by or under contract to, the laboratory. Contracted personnel, when used, must meet competency standards of the laboratory and work in accordance to the laboratory's quality system.



### 17.2 Education and Experience Requirements for Technical Personnel

The laboratory makes every effort to hire analytical staffs that possess a college degree (AA, BA, BS) in an applied science with some chemistry in the curriculum. Exceptions can be made based upon the individual's experience and ability to learn. Selection of qualified candidates for laboratory employment begins with documentation of minimum education, training, and experience prerequisites needed to perform the prescribed task. Minimum education and training requirements for TestAmerica employees are outlined in job descriptions and are generally summarized for analytical staff in the table below.

The laboratory maintains job descriptions for all personnel who manage, perform or verify work affecting the quality of the environmental testing the laboratory performs. Job Descriptions are located on the TestAmerica intranet site's Human Resources web-page (Also see Section 4 for position descriptions/responsibilities).

Experience and specialized training are occasionally accepted in lieu of a college degree (basic lab skills such as using a balance, colony counting, aseptic or quantitation techniques, etc., are also considered).

As a general rule for analytical staff:

Specialty	Education	Experience
Extractions, Digestions, some electrode methods (pH, DO, Redox, etc.), or Titrimetric and Gravimetric Analyses	H.S. Diploma	On the job training (OJT)
GFAA, CVAA, FLAA, Single component or short list Chromatography (e.g., Fuels, BTEX-GC, IC	A college degree in an applied science or 2 years of college and at least 1 year of college chemistry	Or 2 years prior analytical experience is required
ICP, ICPMS, Long List or complex chromatography (e.g., Pesticides, PCB, Herbicides, HPLC, etc.), GCMS	A college degree in an applied science or 2 years of college chemistry	or 5 years of prior analytical experience
Spectra Interpretation	A college degree in an applied science or 2 years of college chemistry	And 2 years relevant experience Or 5 years of prior analytical experience
Technical Managers – <b>General</b>	Bachelors Degree in an applied science or engineering with 24 semester hours in chemistry  An advanced (MS, PhD.) degree may substitute for one year of experience	And 2 years experience in environmental analysis of representative analytes for which they will oversee

Specialty	Education	Experience
Technical Managers – <b>Wet Chem</b> only (no advanced instrumentation)	Associates degree in an applied science or engineering or 2 years of college with 16 semester hours in chemistry	And 2 years relevant experience

When an analyst does not meet these requirements, they can perform a task under the direct supervision of a qualified analyst, peer reviewer or Technical Manager, and are considered an analyst in training. The person supervising an analyst in training is accountable for the quality of the analytical data and must review and approve data and associated corrective actions.

### 17.3 Training

The laboratory is committed to furthering the professional and technical development of employees at all levels.

Orientation to the laboratory's policies and procedures, in-house method training, and employee attendance at outside training courses and conferences all contribute toward employee proficiency. Below are examples of various areas of required employee training:

Required Training	Time Frame	Employee Type
Environmental Health & Safety	Prior to lab work	All
Ethics – New Hires	1 week of hire	All
Ethics – Comprehensive	90 days of hire	All
Data Integrity	30 days of hire	Technical and PMs
Quality Assurance	90 days of hire	All
Ethics – Comprehensive Refresher	Annually	All
Initial Demonstration of Capability (DOC)	Prior to unsupervised method performance	Technical

The laboratory maintains records of relevant authorization/competence, education, professional qualifications, training, skills and experience of technical personnel (including contracted personnel) as well as the date that approval/authorization was given. These records are kept on file at the laboratory. Also refer to "Demonstration of Capability" in Section 19.

The training of technical staff is kept up to date by:

- Each employee must have documentation in their training file that they have read, understood and agreed to follow the most recent version of the laboratory QA Manual and SOPs in their area of responsibility. This documentation is updated as SOPs are updated.
- Documentation from any training courses or workshops on specific equipment, analytical techniques or other relevant topics are maintained in their training file.
- Documentation of proficiency (refer to Section 19).

- An Ethics Agreement signed by each staff member (renewed each year) and evidence of annual ethics training.
- A Confidentiality Agreement signed by each staff member signed at the time of employment.
- Human Resources maintains documentation and attestation forms on employment status & records; benefit programs; timekeeping/payroll; and employee conduct (e.g., ethics violations). This information is maintained in the employee's secured personnel file.

Evidence of successful training could include such items as:

- Adequate documentation of training within operational areas, including one-on-one technical training for individual technologies, and particularly for people cross-trained.
- Analysts knowledge to refer to QA Manual for quality issues.
- Analysts following SOPs, i.e., practice matches SOPs.
- Analysts regularly communicate to supervisors and QA if SOPs need revision, rather than waiting for auditors to find problems.

Further details of the laboratory's training program are described in the Employee Training SOP (TA-QA-0608).

#### 17.4 Data Integrity and Ethics Training Program

Establishing and maintaining a high ethical standard is an important element of a Quality System. Ethics and data integrity training is integral to the success of TestAmerica and is provided for each employee at TestAmerica. It is a formal part of the initial employee orientation within 1 week of hire followed by technical data integrity training within 30 days, comprehensive training within 90 days, and an annual refresher for all employees. Senior management at each facility performs the ethics training for their staff.

In order to ensure that all personnel understand the importance TestAmerica places on maintaining high ethical standards at all times; TestAmerica has established a Corporate Ethics Policy (Policy No. CW-L-P-004) and an Ethics Statement. All initial and annual training is documented by signature on the signed Ethics Statement demonstrating that the employee has participated in the training and understands their obligations related to ethical behavior and data integrity.

Violations of this Ethics Policy will not be tolerated. Employees who violate this policy will be subject to disciplinary actions up to and including termination. Criminal violations may also be referred to the Government for prosecution. In addition, such actions could jeopardize TestAmerica's ability to do work on Government contracts, and for that reason, TestAmerica has a Zero Tolerance approach to such violations.

Employees are trained as to the legal and environmental repercussions that result from data misrepresentation. Key topics covered in the presentation include:

- Organizational mission and its relationship to the critical need for honesty and full disclosure in all analytical reporting.
- Ethics Policy

- How and when to report ethical/data integrity issues. Confidential reporting.
- Record keeping.
- Discussion regarding data integrity procedures.
- Specific examples of breaches of ethical behavior (e.g. peak shaving, altering data or computer clocks, improper macros, etc., accepting/offering kickbacks, illegal accounting practices, unfair competition/collusion)
- Internal monitoring. Investigations and data recalls.
- Consequences for infractions including potential for immediate termination, debarment, or criminal prosecution.
- Importance of proper written narration / data qualification by the analyst and project manager with respect to those cases where the data may still be usable but are in one sense or another partially deficient.

Additionally, a data integrity hotline (1-800-736-9407) is maintained by TestAmerica and administered by the Corporate Quality Department.

## SECTION 18. ACCOMMODATIONS AND ENVIRONMENTAL CONDITIONS

### 18.1 Overview

The *fixed* laboratory is a 20,000 ft<sup>2</sup> secure laboratory facility with controlled access and designed to accommodate an efficient workflow and to provide a safe and comfortable work environment for employees. All visitors sign in and are escorted by laboratory personnel. Access is controlled by various measures.

*The mobile laboratory is a 23 ft van or client provided work space designed to accommodate an efficient workflow and to provide a safe and comfortable work environment for employees.*

The *fixed and mobile* laboratories are equipped with structural safety features. Each employee is familiar with the location, use, and capabilities of general and specialized safety features associated with their workplace. The laboratory provides and requires the use of protective equipment including safety glasses, protective clothing, gloves, etc., OSHA and other regulatory agency guidelines regarding required amounts of bench and fume hood space, lighting, ventilation (temperature and humidity controlled at *fixed laboratory*), access, and safety equipment are met or exceeded.

Traffic flow through sample preparation and analysis areas is minimized to reduce the likelihood of contamination. Adequate floor space and bench top area is provided to allow unencumbered sample preparation and analysis space. Sufficient space is also provided for storage of reagents and media, glassware, and portable equipment. Ample space is also provided for refrigerated sample storage before analysis and archival storage of samples after analysis. The *fixed* laboratory HVAC and deionized water systems are designed to minimize potential trace contaminants.

The *fixed* laboratory is separated into specific areas for sample receiving, sample preparation, volatile organic sample analysis, non-volatile organic sample analysis, inorganic sample analysis and administrative functions.

## 18.2 Environment

Laboratory accommodation, test areas, energy sources, lighting are adequate to facilitate proper performance of tests. The *fixed* facility is equipped with heating, ventilation, and air conditioning (HVAC) systems appropriate to the needs of environmental testing performed at this laboratory.

The environment in which these activities are undertaken does not invalidate the results or adversely affect the required accuracy of any measurements.

The laboratory provides for the effective monitoring, control and recording of environmental conditions that may affect the results of environmental tests as required by the relevant specifications, methods, and procedures. Such environmental conditions include temperature levels in the laboratory.

When any of the method or regulatory required environmental conditions change to a point where they may adversely affect test results, analytical testing will be discontinued until the environmental conditions are returned to the required levels.

Environmental conditions of the facility housing the computer network and LIMS are regulated to protect against raw data loss.

## 18.3 Work Areas

There is effective separation between neighboring areas when the activities therein are incompatible with each other. Examples include:

- Volatile organic chemical handling areas, including sample preparation and waste disposal, and volatile organic chemical analysis areas.

Access to and use of all areas affecting the quality of analytical testing is defined and controlled by secure access to the laboratory building as described below in the Building Security section.

Adequate measures are taken to ensure good housekeeping in the laboratory and to ensure that any contamination does not adversely affect data quality. These measures include regular cleaning to control dirt and dust within the laboratory. Work areas are available to ensure an unencumbered work area. Work areas include:

- Access and entryways to the laboratory.
- Sample receipt areas.
- Sample storage areas.
- Chemical and waste storage areas.
- Data handling and storage areas.

- Sample processing areas.
- Sample analysis areas.

## 18.4 Floor Plan

The *fixed* laboratory floor plan can be found in Appendix 1.

## 18.5 Building Security

Building/Van keys are distributed to employees as necessary.

Visitors to the laboratory sign in and out in a visitor's logbook. A visitor is defined as any person who visits the laboratory who is not an employee of the laboratory. In addition to signing into the laboratory, the Environmental, Health and Safety Manual contains requirements for visitors and vendors. There are specific safety forms that must be reviewed and signed. Visitors (with the exception of company employees) are escorted by laboratory personnel at all times, or the location of the visitor is noted in the visitor's logbook.

## SECTION 19. TEST METHODS AND METHOD VALIDATION

### 19.1 Overview

The laboratory uses methods that are appropriate to meet our clients' requirements and that are within the scope of the laboratory's capabilities. These include sampling, handling, transport, storage and preparation of samples, and, where appropriate, an estimation of the measurement of uncertainty as well as statistical techniques for analysis of environmental data.

Instructions are available in the laboratory for the operation of equipment as well as for the handling and preparation of samples. All instructions, Standard Operating Procedures (SOPs), reference methods and manuals relevant to the working of the laboratory are readily available to all staff. Deviations from published methods are documented (with justification) in the laboratory's approved SOPs. SOPs are submitted to clients for review at their request. Significant deviations from published methods require client approval and regulatory approval where applicable.

### 19.2 Standard Operating Procedures (SOPS)

The laboratory maintains SOPs that accurately reflect all phases of the laboratory such as assessing data integrity, corrective actions, handling customer complaints as well as all analytical methods and sampling procedures. The method SOPs are derived from the most recently promulgated/approved, published methods and are specifically adapted to the laboratory facility. Modifications or clarifications to published methods are clearly noted in the SOPs. All SOPs are controlled in the laboratory.

- All SOPs contain a revision number, effective date, and appropriate approval signatures. Controlled copies are available to all staff.

- Procedures for writing an SOP are incorporated by reference to the laboratory's SOP TA-QA-0500 Standard Operating Procedures.
- SOPs are reviewed at a minimum of every 2 years (annually for Drinking Water and DoD SOPs), and where necessary, revised to ensure continuing suitability and compliance with applicable requirements.

### 19.3 Laboratory Methods Manual

For each test method, the laboratory shall have available the published referenced method as well as the laboratory developed SOP.

**Note:** If more stringent standards or requirements are included in a mandated test method or regulation than those specified in this manual, the laboratory shall demonstrate that such requirements are met. If it is not clear which requirements are more stringent, the standard from the method or regulation is to be followed. Any exceptions or deviations from the referenced methods or regulations are noted in the specific analytical SOP.

The laboratory maintains an SOP Index for both technical and non-technical SOPs. Technical SOPs are maintained to describe a specific test method. Non-technical SOPs are maintained to describe functions and processes not related to a specific test method.

### 19.4 Selection of Methods

Since numerous methods and analytical techniques are available, continued communication between the client and laboratory is imperative to assure the correct methods are utilized. Once client methodology requirements are established, this and other pertinent information is summarized by the Project Manager. These mechanisms ensure that the proper analytical methods are applied when the samples arrive for log-in. For non-routine analytical services (e.g., special matrices, non-routine compound lists), the method of choice is selected based on client needs and available technology. The methods selected should be capable of measuring the specific parameter of interest, in the concentration range of interest, and with the required precision and accuracy.

#### 19.4.1 Sources of Methods

Routine analytical services are performed using standard EPA-approved methodology. In some cases, modification of standard approved methods may be necessary to provide accurate analyses of particularly complex matrices. When the use of specific methods for sample analysis is mandated through project or regulatory requirements, only those methods shall be used.

When clients do not specify the method to be used or methods are not required, the methods used will be clearly validated and documented in an SOP and available to clients and/or the end user of the data.

The analytical methods used by the laboratory are those currently accepted and approved by the U. S. EPA and the state or territory from which the samples were collected. Reference methods include:

- Guidelines Establishing Test Procedures for the Analysis of Pollutants Under the Clean Water Act, and Appendix A-C: 40 CFR Part 136, USEPA Office of Water, Revised as of July 1, 1995, Appendix A to Part 136 - Methods for Organic Chemical Analysis of Municipal and Industrial Wastewater (EPA 600 Series)
- Methods for Chemical Analysis of Water and Wastes, EPA 600 (4-79-020), 1983.
- Methods for the Determination of Inorganic Substances in Environmental Samples, EPA-600/R-93/100, August 1993.
- Methods for the Determination of Metals in Environmental Samples, EPA/600/4-91/010, June 1991. Supplement I: EPA-600/R-94/111, May 1994.
- Methods for the Determination of Organic Compounds in Drinking Water, EPA-600/4-88-039, December 1988, Revised, July 1991, Supplement I, EPA-600-4-90-020, July 1990, Supplement II, EPA-600/R-92-129, August 1992, Supplement III EPA/600/R-95/131 - August 1995 (EPA 500 Series) (EPA 500 Series methods)
- Technical Notes on Drinking Water Methods, EPA-600/R94-173, October 1994
- NIOSH Manual of Analytical Methods, 4<sup>th</sup> ed., August 1994.
- Standard Methods for the Examination of Water and Wastewater, 18<sup>th</sup>/19<sup>th</sup> /20<sup>th</sup>/ on-line edition; Eaton, A.D. Clesceri, L.S. Greenberg, A.E. Eds; American Water Works Association, Water Pollution Control Federation, American Public Health Association: Washington, D.C.
- Test Methods for Evaluating Solid Waste Physical/Chemical Methods (SW846), Third Edition, September 1986, Final Update I, July 1992, Final Update IIA, August 1993, Final Update II, September 1994; Final Update IIB, January 1995; Final Update III, December 1996; Final Update IV, January 2008.
- Annual Book of ASTM Standards, American Society for Testing & Materials (ASTM), Philadelphia, PA.
- National Status and Trends Program, National Oceanographic and Atmospheric Administration, Volume I-IV, 1985-1994.
- Manual for the Certification of Laboratories Analyzing Drinking Water (EPA 815-R-05-004, January 2005).
- Code of Federal Regulations (CFR) 40, Parts 136, 141, 172, 173, 178, 179 and 261
- Plumb, Jr., R.H. 1981, Procedures for Handling and Chemical Analysis of Sediment and Water Samples, Technical Report EPA/C E-81-1. US Army Engineering Waterways Experiment Station, Vicksburg, MS.
- Recommended Protocols for Measuring Conventional Sediment Variables in Puget Sound, April 2003.
- Recommended Guidelines for Measuring Organic Compounds in Puget Sound Waters, Sediment and Tissue Samples, April 1997.
- Recommended Guidelines for Measuring Metals in Puget Sound Waters, Sediment and Tissue Samples, April 1997.

The laboratory reviews updated versions to all the aforementioned references for adaptation based upon capabilities, instrumentation, etc., and implements them as appropriate. As such, the laboratory strives to perform only the latest versions of each approved method as regulations allow or require.

Other reference procedures for non-routine analyses may include methods established by specific states (e.g., Underground Storage Tank methods), ASTM or equipment manufacturers.

## Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 7  
Effective Date: 4/23/2014  
Page 79 of 165

Sample type, source, and the governing regulatory agency requiring the analysis will determine the method utilized.

The laboratory shall inform the client when a method proposed by the client may be inappropriate or out of date. After the client has been informed, and they wish to proceed contrary to the laboratory's recommendation, it will be documented.

### 19.4.2 Demonstration of Capability

Before the laboratory may institute a new method and begin reporting results, the laboratory shall confirm that it can properly operate the method. In general, this demonstration does not test the performance of the method in real world samples, but in an applicable and available clean matrix sample. If the method is for the testing of analytes that are not conducive to spiking, demonstration of capability may be performed on quality control samples.

A demonstration of capability (Analyst DOC, Lab SOP # TA-QA-0617) is performed whenever there is a change in instrument type (e.g., new instrumentation), method or personnel (e.g., analyst hasn't performed the test within the last 12 months).

The initial demonstration of capability must be thoroughly documented and approved by the Technical Manager and QA Manager prior to independently analyzing client samples. All associated documentation must be retained in accordance with the laboratories archiving procedures.

The laboratory must have an approved SOP, demonstrate satisfactory performance, and conduct an MDL study (when applicable). There may be other requirements as stated within the published method or regulations (i.e., retention time window study).

**Note:** In some instances, a situation may arise where a client requests that an unusual analyte be reported using a method where this analyte is not normally reported. If the analyte is being reported for regulatory purposes, the method must meet all procedures outlined within this QA Manual (SOP, MDL, and Demonstration of Capability). If the client states that the information is not for regulatory purposes, the result may be reported as long as the following criteria are met:

- The instrument is calibrated for the analyte to be reported using the criteria for the method and ICV/CCV criteria are met (unless an ICV/CCV is not required by the method or criteria are per project DQOs).
- The laboratory's nominal or default reporting limit (RL) is equal to the quantitation limit (QL), must be at or above the lowest non-zero standard in the calibration curve and must be reliably determined. Project RLs are client specified reporting levels which may be higher than the QL. Results reported below the QL must be qualified as estimated values. Also see Section 19.6.1.3, Relationship of Limit of Detection (LOD) to Quantitation Limit (QL).
- The client request is documented and the lab informs the client of its procedure for working with unusual compounds. The final report must be footnoted: Reporting Limit based on the low standard of the calibration curve.

## Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 7  
Effective Date: 4/23/2014  
Page 80 of 165

### 19.4.3 Initial Demonstration of Capability (IDOC) Procedures

**19.4.3.1** The spiking standard used must be prepared independently from those used in instrument calibration.

**19.4.3.2** The analyte(s) shall be diluted in a volume of clean matrix sufficient to prepare four aliquots at the concentration specified by a method or the laboratory SOP.

**19.4.3.3** At least four aliquots shall be prepared (including any applicable clean-up procedures) and analyzed according to the test method (either concurrently or over a period of days).

**19.4.3.4** Using all of the results, calculate the mean recovery in the appropriate reporting units and the standard deviations for each parameter of interest.

**19.4.3.5** When it is not possible to determine the mean and standard deviations, such as for presence, absence and logarithmic values, the laboratory will assess performance against criteria described in the Method SOP.

**19.4.3.6** Compare the information obtained above to the corresponding acceptance criteria for precision and accuracy in the test method (if applicable) or in laboratory generated acceptance criteria (LCS or interim criteria) if there is no mandatory criteria established. If any one of the parameters do not meet the acceptance criteria, the performance is unacceptable for that parameter.

**19.4.3.7** When one or more of the tested parameters fail at least one of the acceptance criteria, the analyst must proceed according to either option listed below:

- Locate and correct the source of the problem and repeat the test for all parameters of interest beginning with 19.4.3.3 above.
- Beginning with 19.4.3.3 above, repeat the test for all parameters that failed to meet criteria. Repeated failure, however, will confirm a general problem with the measurement system. If this occurs, locate and correct the source of the problem and repeat the test for all compounds of interest beginning with 19.4.3.1 above.

**Note:** Results of successive LCS analyses can be used to fulfill the DOC requirement.

A certification statement (refer to Figure 19-1 as an example) shall be used to document the completion of each initial demonstration of capability. A copy of the certification is archived in the analyst's training folder.

Methods on line prior to the effective date of this Section shall be updated to the procedures outlined above as new analysts perform their demonstration of capability. A copy of the new record will replace that which was used for documentation in the past. At a minimum, the precision and accuracy of four mid-level laboratory control samples must have been compared to the laboratory's quality control acceptance limits.

**19.5 Laboratory Developed Methods and Non-Standard Methods**

Any new method developed by the laboratory must be fully defined in an SOP and validated by qualified personnel with adequate resources to perform the method. Method specifications and the relation to client requirements must be clearly conveyed to the client if the method is a non-standard method (not a published or routinely accepted method). The client must also be in agreement to the use of the non-standard method.

**19.6 Validation of Methods**

Validation is the confirmation by examination and the provision of objective evidence that the particular requirements for a specific intended use are fulfilled.

All non-standard methods, laboratory designed/developed methods, standard methods used outside of their scope, and major modifications to published methods must be validated to confirm they are fit for their intended use. The validation will be as extensive as necessary to meet the needs of the given application. The results are documented with the validation procedure used and contain a statement as to the fitness for use.

**19.6.1 Method Validation and Verification Activities for All New Methods**

While method validation can take various courses, the following activities can be required as part of method validation. Method validation records are designated QC records and are archived accordingly.

**19.6.1.1 Determination of Method Selectivity**

Method selectivity is the demonstrated ability to discriminate the analyte(s) of interest from other compounds in the specific matrix or matrices from other analytes or interference. In some cases to achieve the required selectivity for an analyte, a confirmation analysis is required as part of the method.

**19.6.1.2 Determination of Method Sensitivity**

Sensitivity can be both estimated and demonstrated. Whether a study is required to estimate sensitivity depends on the level of method development required when applying a particular measurement system to a specific set of samples. Where estimations and/or demonstrations of sensitivity are required by regulation or client agreement, such as the procedure in 40 CFR Part 136 Appendix B, under the Clean Water Act, these shall be followed.

**19.6.1.3 Relationship of Limit of Detection (LOD) to the Quantitation Limit (QL)**

An important characteristic of expression of sensitivity is the difference in the LOD and the QL. The LOD is the minimum level at which the presence of an analyte can be reliably concluded. The QL is the minimum concentration of analyte that can be quantitatively determined with acceptable precision and bias. For most instrumental measurement systems, there is a region where semi-quantitative data is generated around the LOD (both above and below the estimated MDL or LOD) and below the QL. In this region, detection of an analyte may be confirmed but quantification of the analyte is unreliable within the accuracy and precision guidelines of the measurement system. When an analyte is detected below the QL, and the

presence of the analyte is confirmed by meeting the qualitative identification criteria for the analyte, the analyte can be reliably reported, but the amount of the analyte can only be estimated. If data is to be reported in this region, it must be done so with a qualification that denotes the semi-quantitative nature of the result.

**19.6.1.4 Determination of Interferences**

A determination that the method is free from interferences in a blank matrix is performed.

**19.6.1.5 Determination of Range**

Where appropriate to the method, the quantitation range is determined by comparison of the response of an analyte in a curve to established or targeted criteria. Generally the upper quantitation limit is defined by highest acceptable calibration concentration. The lower quantitation limit or QL cannot be lower than the lowest non-zero calibration level, and can be constrained by required levels of bias and precision.

**19.6.1.6 Determination of Accuracy and Precision**

Accuracy and precision studies are generally performed using replicate analyses, with a resulting percent recovery and measure of reproducibility (standard deviation, relative standard deviation) calculated and measured against a set of target criteria.

**19.6.1.7 Documentation of Method**

The method is formally documented in an SOP. If the method is a minor modification of a standard laboratory method that is already documented in an SOP, an SOP Attachment describing the specific differences in the new method is acceptable in place of a separate SOP.

**19.6.1.8 Continued Demonstration of Method Performance**

Continued demonstration of Method Performance is addressed in the SOP. Continued demonstration of method performance is generally accomplished by batch specific QC samples such as LCS, method blanks or PT samples.

**19.7 Method Detection Limits (MDL) / Limits of Detection (LOD)**

Method detection limits (MDL) are initially determined in accordance with 40 CFR Part 136, Appendix B or alternatively by other technically acceptable practices that have been accepted by regulators. MDL is also sometimes referred to as Limit of Detection (LOD). The MDL theoretically represents the concentration level for each analyte within a method at which the Analyst is 99% confident that the true value is not zero. The MDL is determined for each analyte initially during the method validation process and updated as required in the analytical methods, whenever there is a significant change in the procedure or equipment, or based on project specific requirements. Generally, the analyst prepares at least seven replicates of solution spiked at one to five times the estimated method detection limit (most often at the lowest standard in the calibration curve) into the applicable matrix with all the analytes of interest. Each of these aliquots is extracted (including any applicable clean-up procedures) and analyzed in the same manner as the samples. Where possible, the seven replicates should be analyzed over 2-4 days to provide



a more realistic MDL. To allow for some flexibility, this low level standard may be analyzed every batch or every week or some other frequency rather than doing the study all at once. In addition, a larger number of data points may be used if the appropriate t-value multiplier is used.

Refer to the Corporate SOP No. CA-Q-S-006 or the laboratory's SOP No. TA-QA-0602 for details on the laboratory's MDL process.

#### 19.8 Instrument Detection Limits (IDL)

The IDL is sometimes used to assess the reasonableness of the MDLs or in some cases required by the analytical method or program requirements. IDLs are most used in metals analyses, but may be useful in demonstration of instrument performance in other areas.

IDLs are calculated to determine an instrument's sensitivity independent of any preparation method. IDLs are calculated either using 7 replicate spike analyses, like MDL but without sample preparation, or by the analysis of 10 instrument blanks and calculating 3 x the absolute value of the standard deviation.

If IDL is > than the MDL, it may be used as the reported MDL.

#### 19.9 Verification of Detection and Reporting Limits

Once the MDL is determined, it must be verified on each instrument used for the given method. TestAmerica defines the DoD QSM Detection Limit (DL) as being equal to the MDL. TestAmerica also defines the DoD QSM Limit of Detection (LOD) as being equal to the lowest concentration standard that successfully verifies the MDL, also referred to as the MDLV standard. MDL and MDLV standards are extracted/digested and analyzed through the entire analytical process. The MDL and MDLV determinations do not apply to methods that are not readily spiked (e.g. pH, turbidity, etc.) or where the lab does not report to the MDL. If the MDLV standard is not successful, then the laboratory will redevelop their MDL or perform and pass two consecutive MDLVs at a higher concentration and set the LOD at the higher concentration. Initial and quarterly verification is required for all methods listed in the laboratory's DoD ELAP Scope of Accreditation. Refer to the laboratory SOP TA-QA-0602 Method Detection Limits (MDLs/DLs) for further details.

The laboratory quantitation limit is equivalent to the DoD Limit of Quantitation (LOQ), which is at a concentration equal to or greater than the lowest non-zero calibration standard. The DoD QSM requires the laboratory to perform an initial characterization of the bias and precision at the LOQ and quarterly LOQ verifications thereafter. If the quarterly verification results are not consistent with three-standard deviation confidence limits established initially, then the bias and precision will be reevaluated and clients contacted for any on-going projects. For DoD projects, TestAmerica makes a distinction between the Reporting Limit (RL) and the LOQ. The RL is a level at or above the LOQ that is used for specific project reporting purposes, as agreed to between the laboratory and the client. The RL cannot be lower than the LOQ concentration, but may be higher.

#### 19.10 Retention Time Windows

Most organic analyses and some inorganic analyses use chromatography techniques for qualitative and quantitative determinations. For every chromatography analysis or as specific in the reference method, each analyte will have a specific time of elution from the column to the detector. This is known as the analyte's retention time. The variance in the expected time of elution is defined as the retention time window. As the key to analyte identification in chromatography, retention time windows must be established on every column for every analyte used for that method. These records are kept with the files associated with an instrument for later quantitation of the analytes. Complete details are available in the laboratory SOPs.

#### 19.11 Evaluation of Selectivity

The laboratory evaluates selectivity by following the checks within the applicable analytical methods, which include mass spectral tuning, second column confirmation, ICP interelement interference checks, chromatography retention time windows, sample blanks, spectrochemical, atomic absorption profiles, co-precipitation evaluations and specific electrode response factors.

#### 19.12 Estimation of Uncertainty of Measurement

**19.12.1** Uncertainty is "a parameter associated with the result of a measurement, that characterizes the dispersion of the values that could reasonably be attributed to the measurand" (as defined by the International Vocabulary of Basic and General Terms in Metrology, ISO Geneva, 1993, ISBN 92-67-10175-1). Knowledge of the uncertainty of a measurement provides additional confidence in a result's validity. Its value accounts for all the factors which could possibly affect the result, such as adequacy of analyte definition, sampling, matrix effects and interferences, climatic conditions, variances in weights, volumes, and standards, analytical procedure, and random variation. Some national accreditation organizations require the use of an "expanded uncertainty": the range within which the value of the measurand is believed to lie within at least a 95% confidence level with the coverage factor k=2.

**19.12.2** Uncertainty is not error. Error is a single value, the difference between the true result and the measured result. On environmental samples, the true result is never known. The measurement is the sum of the unknown true value and the unknown error. Unknown error is a combination of systematic error, or bias, and random error. Bias varies predictably, constantly, and independently from the number of measurements. Random error is unpredictable, assumed to be Gaussian in distribution, and reducible by increasing the number of measurements.

**19.12.3** The minimum uncertainty associated with results generated by the laboratory can be determined by using the Laboratory Control Sample (LCS) accuracy range for a given analyte. The LCS limits are used to assess the performance of the measurement system since they take into consideration all of the laboratory variables associated with a given test over time (except for variability associated with the sampling and the variability due to matrix effects). The percent recovery of the LCS is compared either to the method-required LCS accuracy limits or to the statistical, historical, in-house LCS accuracy limits.

**19.12.4** To calculate the uncertainty for the specific result reported, multiply the result by the decimal of the lower end of the LCS range percent value for the lower end of the uncertainty



range, and multiply the result by the decimal of the upper end of the LCS range percent value for the upper end of the uncertainty range. These calculated values represent uncertainties at approximately the 99% confidence level with a coverage factor of  $k = 3$ . As an example, for a reported result of 1.0 mg/l with an LCS recovery range of 50 to 150%, the estimated uncertainty in the result would be 1.0 +/- 0.5 mg/l.

**19.12.5** In the case where a well recognized test method specifies limits to the values of major sources of uncertainty of measurement (e.g., 524.2, 525, etc.) and specifies the form of presentation of calculated results, no further discussion of uncertainty is required.

#### **19.13 Sample Reanalysis Guidelines**

Because there is a certain level of uncertainty with any analytical measurement, a sample re-preparation (where appropriate) and subsequent analysis (hereafter referred to as 'reanalysis') may result in either a higher or lower value from an initial sample analysis. There are also variables that may be present (e.g., sample homogeneity, analyte precipitation over time, etc.) that may affect the results of a reanalysis. Based on the above comments, the laboratory will reanalyze samples at a client's request with the following caveats. Client specific Contractual Terms & Conditions for reanalysis protocols may supersede the following items.

- Homogenous samples: If a reanalysis agrees with the original result to within the RPD limits for MS/MSD or Duplicate analyses, or within  $\pm 1$  reporting limit for samples  $\leq 5x$  the reporting limit, the original analysis will be reported. At the client's request, both results may be reported on the same report but not on two separate reports.
- If the reanalysis does not agree (as defined above) with the original result, then the laboratory will investigate the discrepancy and reanalyze the sample a third time for confirmation if sufficient sample is available.
- Any potential charges related to reanalysis are discussed in the contract terms and conditions or discussed at the time of the request. The client will typically be charged for reanalysis unless it is determined that the lab was in error.
- Due to the potential for increased variability, reanalysis may not be applicable to Non-homogenous, Encore, and Sodium Bisulfate preserved samples. See the Area Technical Manager or Laboratory Director if unsure.

#### **19.14 Control of Data**

The laboratory has policies and procedures in place to ensure the authenticity, integrity, and accuracy of the analytical data generated by the laboratory.

##### **19.14.1 Computer and Electronic Data Related Requirements**

The three basic objectives of our computer security procedures and policies are shown below. More detail is outlined in SOP TA-IT-0523 Data Integrity and Security. The laboratory is currently running the TALS system which is a custom in-house developed LIMS system that has been highly customized to meet the needs of the laboratory. It is referred to as LIMS for the

remainder of this section. The LIMS utilizes Sequel Server which is an industry standard relational database platform. It is referred to as Database for the remainder of this section.

**19.14.1.1 Maintain the Database Integrity:** Assurance that data is reliable and accurate through data verification (review) procedures, password-protecting access, anti-virus protection, data change requirements, as well as an internal LIMS permissions procedure.

- LIMS Database Integrity is achieved through data input validation, internal user controls, and data change requirements.
- Spreadsheets and other software developed in-house must be verified with documentation through hand calculations prior to use. Cells containing calculations must be lock-protected and controlled.
- Instrument hardware and software adjustments are safeguarded through maintenance logs, audit trails and controlled access.

**19.14.1.2 Ensure Information Availability:** Protection against loss of information or service is ensured through scheduled back-ups, stable file server network architecture, secure storage of media, line filter, Uninterruptible Power Supply (UPS), and maintaining older versions of software as revisions are implemented.

**19.14.1.3 Maintain Confidentiality:** Ensure data confidentiality through physical access controls such as password protection or website access approval when electronically transmitting data.

##### **19.14.2 Data Reduction**

The complexity of the data reduction depends on the analytical method and the number of discrete operations involved (e.g., extractions, dilutions, instrument readings and concentrations). The analyst calculates the final results from the raw data or uses appropriate computer programs to assist in the calculation of final reportable values.

For manual data entry, e.g., Wet Chemistry, the data is reduced by the analyst and then verified by the Technical Manager or alternate analyst prior to updating the data in LIMS. The spreadsheets, or any other type of applicable documents, are signed by both the analyst and alternate reviewer to confirm the accuracy of the manual entry(s).

Manual integration of peaks will be documented and reviewed and the raw data will be flagged in accordance with the TestAmerica Corporate SOP No. CA-Q-S-002, Acceptable Manual Integration Practices.

Analytical results are reduced to appropriate concentration units specified by the analytical method, taking into account factors such as dilution, sample weight or volume, etc. Blank correction will be applied only when required by the method or per manufacturer's indication; otherwise, it should not be performed. Calculations are independently verified by appropriate laboratory staff. Calculations and data reduction steps for various methods are summarized in the respective analytical SOPs or program requirements.

**19.14.2.1** All raw data must be retained in the worklist folder, computer file (if appropriate), and/or runlog. All criteria pertinent to the method must be recorded. The

documentation is recorded at the time observations or calculations are made and must be signed or initialed/dated (month/day/year). It must be easily identifiable who performed which tasks if multiple people were involved.

- 19.14.2.2** In general, concentration results are reported in milligrams per liter (mg/L) or micrograms per liter (µg/l) for liquids and milligrams per kilogram (mg/Kg) or micrograms per kilogram (µg/Kg) for solids. For values greater than 10,000 mg/L, results can be reported in percent, i.e., 10,000 mg/l = 1%. Units are defined in each lab SOP.
- 19.14.2.3** In reporting, the analyst or the instrument output records the raw data result using values of known certainty plus one uncertain digit. If final calculations are performed external to LIMS, the results should be entered in LIMS with at least three significant figures. In general, results are reported to 2 significant figures on the final report.
- 19.14.2.4** For those methods that do not have an instrument printout or an instrumental output compatible with the LIMS System, the raw results and dilution factors are entered directly into LIMS by the analyst, and the software calculates the final result for the analytical report. LIMS has a defined significant figure criterion for each analyte.
- 19.14.2.5** The laboratory strives to import data directly from instruments or calculation spreadsheets to ensure that the reported data are free from transcription and calculation errors. For those analyses with an instrumental output compatible with the LIMS, the raw results and dilution factors are transferred into LIMS electronically after reviewing the quantitation report, and removing unrequested or poor spectrally-matched compounds. The analyst prints a copy of what has been entered to check for errors. This printout and the instrument's printout of calibrations, concentrations, retention times, chromatograms, and mass spectra, if applicable, are retained with the data file. The data file is stored in a monthly folder on the instrument computer; periodically, this file is transferred to the server and, eventually, to a tape file.

**19.14.3** Logbook / Worksheet Use Guidelines

Logbooks and worksheets are filled out 'real time' and have enough information on them to trace the events of the applicable analysis/task. (e.g. calibrations, standards, analyst, sample ID, date, time on short holding time tests, temperatures when applicable, calculations are traceable, etc.)

- Corrections are made following the procedures outlined in Section 12.
- Logbooks are controlled by the QA department. A record is maintained of all logbooks in the lab.
- Unused portions of pages must be "Z"ed out, signed and dated.
- Worksheets are created with the approval of the Technical Director and QA Manager at the facility. The QA Manager controls all worksheets following the procedures in Section 6.

**19.14.4** Review / Verification Procedures

Review procedures are outlined in several SOPs and SOP TA-QA-0635 Procedure for Peer Review of Analytical Data and Project Management Review of Reports to ensure that reported

data are free from calculation and transcription errors and that QC parameters have been reviewed and evaluated before data is reported. The laboratory also utilizes corporate SOP CA-Q-S-002 Acceptable Manual Integration Practices to ensure the authenticity of the data. The general review concepts are discussed below, more specific information can be found in the SOPs.

- 19.14.4.1** The data review process at the laboratory starts at the Sample Control level. Sample Control personnel review chain-of-custody forms and input the sample information and required analyses into a computer LIMS. The Sample Control Supervisor reviews the transaction of the chain-of-custody forms and the inputted information. The Project Managers perform final review of the chain-of-custody forms and inputted information.

- 19.14.4.2** The next level of data review occurs with the Analysts. As results are generated, analysts review their work to ensure that the results generated meet QC requirements and relevant EPA methodologies. The Analysts transfer the data into the LIMS and add data qualifiers if applicable. To ensure data compliance, a different analyst performs a second level of review. Second level review is accomplished by checking reported results against raw data and evaluating the results for accuracy. During the second level review, blank runs, QA/QC check results, initial and continuing calibration results, laboratory control samples, sample data, qualifiers and spike information are evaluated. Where calibration is not required on a daily basis, secondary review of the initial calibration results may be conducted at the time of calibration. Approximately 15% of all sample data from manual methods and from automated methods, all GC/MS spectra and all manual integrations are reviewed. Manual integrations are also electronically reviewed utilizing auditing software to help ensure compliance to ethics and manual integration policies. Issues that deem further review include the following:

- QC data are outside the specified control limits for accuracy and precision
- Reviewed sample data does not match with reported results
- Unusual detection limit changes are observed
- Samples having unusually high results
- Samples exceeding a known regulatory limit
- Raw data indicating some type of contamination or poor technique
- Inconsistent peak integration
- Transcription errors
- Results outside of calibration range

- 19.14.4.3** Unacceptable analytical results may require reanalysis of the samples. Any problems are brought to the attention of the Laboratory Director, Project Manager, Quality Assurance Manager, Technical Manager, or Supervisor for further investigation. Corrective action is initiated whenever necessary.

- 19.14.4.4** The results are then entered or directly transferred into the computer database and a hard copy (or .pdf) is printed for the client.

- 19.14.4.5** As a final review prior to the release of the report, the Project Manager reviews the results for appropriateness and completeness. This review and approval ensures that client requirements have been met and that the final report has been properly completed. The process includes, but is not limited to, verifying that chemical relationships are evaluated, COC is followed, cover letters/ narratives are present, flags are appropriate, and project specific requirements are met.
- 19.14.4.6** Any project that requires a data package is subject to a tertiary data review for transcription errors and acceptable quality control requirements. The Project Manager then signs the final report. The accounting personnel also check the report for any clerical or invoicing errors. When complete, the report is sent out to the client.
- 19.14.4.7** A visual summary of the flow of samples and information through the laboratory, as well as data review and validation, is presented in Figure 19-2.

**19.14.5      Manual Integrations**

Computerized data systems provide the analyst with the ability to re-integrate raw instrument data in order to optimize the interpretation of the data. Though manual integration of data is an invaluable tool for resolving variations in instrument performance and some sample matrix problems, when used improperly, this technique would make unacceptable data appear to meet quality control acceptance limits. Improper re-integrations lead to legally indefensible data, a poor reputation, or possible laboratory decertification. Because guidelines for re-integration of data are not provided in the methods and most methods were written prior to widespread implementation of computerized data systems, the laboratory trains all analytical staff on proper manual integration techniques using TestAmerica's Corporate SOP (CA-Q-S-002).

- 19.14.5.1** The analyst must adjust baseline or the area of a peak in some situations, for example when two compounds are not adequately resolved or when a peak shoulder needs to be separated from the peak of interest. The analyst must use professional judgment and common sense to determine when manual integrating is required. Analysts are encouraged to ask for assistance from a senior analyst or manager when in doubt.
- 19.14.5.2** Analysts shall not increase or decrease peak areas for the sole purpose of achieving acceptable QC recoveries that would have otherwise been unacceptable. The intentional recording or reporting of incorrect information (or the intentional omission of correct information) is against company principals and policy and is grounds for immediate termination.
- 19.14.5.3** Client samples, performance evaluation samples, and quality control samples are all treated equally when determining whether or not a peak area or baseline should be manually adjusted.
- 19.14.5.4** All manual integrations receive a second level review. Manual integrations must be indicated on an expanded scale "after" chromatograms such that the integration performed can be easily evaluated during data review. Expanded scale "before" chromatograms are also required for all manual integrations on QC parameters (calibrations, calibration verifications, laboratory control samples, internal standards,

surrogates, etc.) unless the laboratory has another documented corporate approved procedure in place that can demonstrate an active process for detection and deterrence of improper integration practices.



## SECTION 20. EQUIPMENT and CALIBRATIONS

### 20.1 Overview

The laboratory purchases the most technically advanced analytical instrumentation for sample analyses. Instrumentation is purchased on the basis of accuracy, dependability, efficiency and sensitivity. Each laboratory is furnished with all items of sampling, preparation, analytical testing and measurement equipment necessary to correctly perform the tests for which the laboratory has capabilities. Each piece of equipment is capable of achieving the required accuracy and complies with specifications relevant to the method being performed. Before being placed into use, the equipment (including sampling equipment) is calibrated and checked to establish that it meets its intended specification. The calibration routines for analytical instruments establish the range of quantitation. Calibration procedures are specified in laboratory SOPs. A list of laboratory instrumentation is presented in Table 20-1. A list of software is provided in Table 20-2.

Equipment is only operated by authorized and trained personnel. Manufacturer instructions for equipment use are readily accessible to all appropriate laboratory personnel.

### 20.2 Preventive Maintenance

The laboratory follows a well-defined maintenance program to ensure proper equipment operation and to prevent the failure of laboratory equipment or instrumentation during use. This program of preventive maintenance helps to avoid delays due to instrument failure.

Routine preventive maintenance procedures and frequency, such as cleaning and replacements, should be performed according to the procedures outlined in the manufacturer's manual. Qualified personnel must also perform maintenance when there is evidence of degradation of peak resolution, a shift in the calibration curve, loss of sensitivity, or failure to continually meet one of the quality control criteria.

Table 20-3 lists examples of scheduled routine maintenance. It is the responsibility of each Technical Manager to ensure that instrument maintenance logs are kept for all equipment in his/her department. Preventative maintenance procedures are also outlined in analytical SOPs or instrument manuals. (Note: for some equipment, the log used to monitor performance is also the maintenance log. Multiple pieces of equipment may share the same log as long as it is clear as to which instrument is associated with an entry.)

Instrument maintenance logs are controlled and are used to document instrument problems, instrument repair and maintenance activities. Maintenance logs shall be kept for all major pieces of equipment. Instrument maintenance logs may also be used to specify instrument parameters.

- Documentation must include all major maintenance activities such as contracted preventive maintenance and service and in-house activities such as the replacement of electrical components, lamps, tubing, valves, columns, detectors, cleaning and adjustments.
- Each entry in the instrument log includes the Analyst's initials, the date, a detailed description of the problem (or maintenance needed/scheduled), a detailed explanation of the solution or maintenance performed, and a verification that the equipment is functioning properly (state what was used to determine a return to control. e.g. CCV run on 'date' was acceptable, or

instrument recalibrated on 'date' with acceptable verification, etc.) must also be documented in the instrument records.

- When maintenance or repair is performed by an outside agency, service receipts detailing the service performed can be affixed into the logbooks adjacent to pages describing the maintenance performed. This stapled in page must be signed across the page entered and the logbook so that it is clear that a page is missing if only half a signature is found in the logbook.

If an instrument requires repair (subjected to overloading or mishandling, gives suspect results, or otherwise has shown to be defective or outside of specified limits) it shall be taken out of operation and tagged as out-of-service or otherwise isolated until such a time as the repairs have been made and the instrument can be demonstrated as operational by calibration and/or verification or other test to demonstrate acceptable performance. The laboratory shall examine the effect of this defect on previous analyses.

In the event of equipment malfunction that cannot be resolved, service shall be obtained from the instrument vendor manufacturer, or qualified service technician, if such a service can be tendered. If on-site service is unavailable, arrangements shall be made to have the instrument shipped back to the manufacturer for repair. Back up instruments, which have been approved, for the analysis shall perform the analysis normally carried out by the malfunctioning instrument. If the back up is not available and the analysis cannot be carried out within the needed timeframe, the samples shall be subcontracted.

If an instrument is sent out for service or transferred to another facility, it must be recalibrated and verified (including new initial MDL study) prior to return to lab operations.

### 20.3 Support Equipment

This section applies to all devices that may not be the actual test instrument, but are necessary to support laboratory operations. These include but are not limited to: balances, ovens, refrigerators, freezers, incubators, water baths, temperature measuring devices, thermal/pressure sample preparation devices and volumetric dispensing devices if quantitative results are dependent on their accuracy, as in standard preparation and dispensing or dilution into a specified volume. All raw data records associated with the support equipment are retained to document instrument performance.

#### 20.3.1 Weights and Balances

The accuracy of the balances used in the laboratory is checked every working day, before use. All balances are placed on stable counter tops.

Each balance is checked prior to initial serviceable use with at least two certified ASTM type 1 weights spanning its range of use (weights that have been calibrated to ASTM type 1 weights may also be used for daily verification). ASTM type 1 weights used only for calibration of other weights (and no other purpose) are inspected for corrosion, damage or nicks at least annually and if no damage is observed, they are calibrated at least every 5 years by an outside calibration laboratory. Any weights (including ASTM Type 1) used for daily balance checks or other purposes are recalibrated/recertified annually to NIST standards (this may be done internally if laboratory maintains "calibration only" ASTM type 1 weights).

All balances are serviced annually by a qualified service representative, who supplies the laboratory with a certificate that identifies traceability of the calibration to the NIST standards.

All of this information is recorded in logs, and the recalibration/recertification certificates are kept on file in the QA office. Refer to SOP TA-QA-0014 Selecting and Using Balances for details.

#### **20.3.2 pH, Conductivity, and Turbidity Meters**

The pH meters used in the laboratory are accurate to  $\pm 0.1$  pH units, and have a scale readability of at least 0.05 pH units. The meters automatically compensate for the temperature, and are calibrated with at least two working range buffer solutions before each use.

Conductivity meters are also calibrated before each use with a known standard to demonstrate the meters do not exceed an error of 1% or one umhos/cm.

Turbidity meters are also calibrated before each use. All of this information is documented in logs.

Consult pH and Conductivity, and Turbidity SOPs for further information.

#### **20.3.3 Thermometers**

All thermometers are calibrated on an annual basis with a NIST-traceable thermometer at temperatures bracketing the range of use. IR thermometers, digital probes and thermocouples are calibrated quarterly.

The mercury NIST thermometer is recalibrated every five years (unless thermometer has been exposed to temperature extremes or apparent separation of internal liquid) by an approved outside service and the provided certificate of traceability is kept on file. The NIST thermometer(s) have increments of 1 degree (0.5 degree or less increments are required for drinking water microbiological laboratories), and have ranges applicable to method and certification requirements. The NIST traceable thermometer is used for no other purpose than to calibrate other thermometers.

All of this information is documented in logbooks. Monitoring method-specific temperatures, including incubators, heating blocks, water baths, and ovens, is documented in method-specific logbooks. More information on this subject can be found in the SOP TA-QA-0024 Use, Calibration, and Maintenance of Laboratory Thermometers.

#### **20.3.4 Refrigerators/Freezer Units, Waterbaths, Ovens and Incubators**

The temperatures of all refrigerator units and freezers used for sample and standard storage are monitored each working day. As a DOD-approved lab, sample storage is monitored 7 days a week.

Ovens, water baths and incubators are monitored on days of use.

All of this equipment has a unique identification number, and is assigned a unique thermometer for monitoring.

Sample storage refrigerator temperatures are kept between  $> 0^{\circ}\text{C}$  and  $\leq 6^{\circ}\text{C}$ .

Specific temperature settings/ranges for other refrigerators, ovens, water baths, and incubators can be found in method specific SOPs.

All of this information is documented in Daily Temperature Logbooks and method-specific logbooks.

#### **20.3.5 Autopipettors, Dilutors, and Syringes**

Mechanical volumetric dispensing devices including burettes (except Class A Glassware) are given unique identification numbers and the delivery volumes are verified gravimetrically, at a minimum, on a quarterly basis.

For those dispensers that are not used for analytical measurements, a label is applied to the device stating that it is not calibrated. Any device not regularly verified can not be used for any quantitative measurements. Refer to SOP TA-QA-0016 Volumetric Verification.

Micro-syringes are purchased from Hamilton Company. Each syringe is traceable to NIST. The laboratory keeps on file an "Accuracy and Precision Statement of Conformance" from Hamilton attesting established accuracy.

#### **20.4 Instrument Calibrations**

Calibration of analytical instrumentation is essential to the production of quality data. Strict calibration procedures are followed for each method. These procedures are designed to determine and document the method detection limits, the working range of the analytical instrumentation and any fluctuations that may occur from day to day.

Sufficient raw data records are retained to allow an outside party to reconstruct all facets of the initial calibration. Records contain, but are not limited to, the following: calibration date, method, instrument, analyst(s) initials or signatures, analysis date, analytes, concentration, response, type of calibration (Avg RF, curve, or other calculations that may be used to reduce instrument responses to concentration.)

Sample results must be quantitated from the initial calibration and may not be quantitated from any continuing instrument calibration verification unless otherwise required by regulation, method or program.

If the initial calibration results are outside of the acceptance criteria, corrective action is performed and any affected samples are reanalyzed if possible. If the reanalysis is not possible, any data associated with an unacceptable initial calibration will be reported with appropriate data qualifiers (refer to Section 12).

**Note:** Instruments are calibrated initially and as needed after that and at least annually.

#### 20.4.1 Calibration Standards

Calibration standards are prepared using the procedures indicated in the Reagents and Standards section of the determinative method SOP. If a reference method does not specify the number of calibration standards, a minimum of 3 calibration points (exception being ICP and ICP/MS methods) will be used.

Standards for instrument calibration are obtained from a variety of sources. All standards are traceable to national or international standards of measurement, or to national or international standard reference materials.

The lowest concentration calibration standard that is analyzed during an initial calibration must be at or below the stated reporting limit for the method based on the final volume of extract (or sample).

The other concentrations define the working range of the instrument/method or correspond to the expected range of concentrations found in actual samples that are also within the working range of the instrument/method. Results of samples not bracketed by initial instrument calibration standards (within calibration range to at least the same number of significant figures used to report the data) must be reported as having less certainty, e.g., defined qualifiers or flags (additional information may be included in the case narrative). The exception to these rules is ICP methods or other methods where the referenced method does not specify two or more standards.

All initial calibrations are verified with a standard obtained from a second source and traceable to a national standard, when available (or vendor certified different lot if a second source is not available). This verification occurs immediately after the calibration curve has been analyzed, and before the analysis of any samples.

##### 20.4.1.1 Calibration Verification

The calibration relationship established during the initial calibration must be verified initially and at least daily as specified in the laboratory method SOPs in accordance with the referenced analytical methods and in the 2009 TNI Standard. The process of calibration verification applies to both external standard and internal standard calibration techniques, as well as to linear and non-linear calibration models. Initial calibration verification is with a standard source secondary (second source standard) to the calibration standards, but continuing calibration verifications may use the same source standards as the calibration curve.

**Note:** The process of calibration verification referred to here is fundamentally different from the approach called "calibration" in some methods. As described in those methods, the calibration factors or response factors calculated during calibration are used to update the calibration factors or response factors used for sample quantitation. This approach, while employed in other EPA programs, amounts to a daily single-point calibration.

All target analytes and surrogates, including those reported as non-detects, must be included in periodic calibration verifications for purposes of retention time confirmation and to demonstrate that calibration verification criteria are being met, i.e., RPD, per 2009 TNI Std. EL-V1M4 Sec. 1.7.2.

All samples must be bracketed by periodic analyses of standards that meet the QC acceptance criteria (e.g., calibration and retention time). The frequency is found in the determinative methods or SOPs.

Generally, the initial calibrations must be verified at the beginning of each 12-hour analytical shift during which samples are analyzed. (Some methods may specify more or less frequent verifications). The 12-hour analytical shift begins with the injection of the calibration verification standard (or the MS tuning standard in MS methods). The shift ends after the completion of the analysis of the last sample, QC, or standard that can be injected within 12 hours of the beginning of the shift.

A continuing instrument calibration verification (CCV) must be repeated at the beginning and, for methods that have quantitation by external calibration models, at the end of each analytical batch. Some methods have more frequent CCV requirements see specific SOPs. Most Inorganic methods require the CCV to be analyzed after every 10 samples or injections, including matrix or batch QC samples.

**Note:** If an internal standard calibration is being used (basically GCMS) then bracketing standards are not required, only daily verifications are needed. The results from these verification standards must meet the calibration verification criteria and the retention time criteria (if applicable).

##### 20.4.1.2 Verification of Linear and Non-Linear Calibrations

Calibration verification for calibrations involves the calculation of the percent drift or the percent difference of the instrument response between the initial calibration and each subsequent analysis of the verification standard. (These calculations are available in the laboratory method SOPs. Verification standards are evaluated based on the % Difference from the average CF or RF of the initial calibration or based on % Drift or % Recovery if a linear or quadratic curve is used.

Regardless of whether a linear or non-linear calibration model is used, if the results of a CCV are outside the established acceptance criteria and analysis of a second consecutive (and immediate) CCV fails to produce results within acceptance criteria, corrective action shall be performed. Once corrective actions have been completed & documented, the laboratory shall demonstrate acceptable instrument / method performance by analyzing two consecutive CCVs, or a new initial instrument calibration shall be performed.

Sample analyses and reporting of data may not occur or continue until the analytical system is calibrated or calibration verified. However, data associated with an unacceptable calibration verification may be fully useable under the following special conditions and reported based upon discussion and approval of the client:

- when the acceptance criteria for the CCV are exceeded high (i.e., high bias) and the associated samples within the batch are non-detects, then those non-detects may be reported with a footnote or case narrative explaining the high bias. Otherwise the samples affected by the unacceptable CCV shall be re-analyzed after a new calibration curve has been established, evaluated and accepted; or



b). when the acceptance criteria for the CCV are exceeded low (i.e., low bias), those sample results may be reported if they exceed a maximum regulatory limit/decision level. Otherwise the samples affected by the unacceptable CCV shall be re-analyzed after a new calibration curve has been established, evaluated and accepted. Alternatively, a reporting limit standard may be analyzed to demonstrate that the laboratory can still support non-detects at their reporting limit.

Samples reported by either of the conditions identified above will be appropriately flagged.

#### 20.5 Tentatively Identified Compounds (TICs) – GC/MS Analysis

For samples containing components not associated with the calibration standards, a library search may be made for the purpose of tentative identification. The necessity to perform this type of identification will be determined by the purpose of the analyses being conducted. Data system library search routines should not use normalization routines that would misrepresent the library or unknown spectra when compared to each other.

**Note:** If the TIC compound is not part of the client target analyte list but is calibrated by the laboratory and is both qualitatively and/or quantitatively identifiable, it should not be reported as a TIC. If the compound is reported on the same form as true TICs, it should be qualified and/or narrated that the reported compound is qualitatively and quantitatively (if verification in control) reported compared to a known standard that is in control (where applicable).

For example, the RCRA permit or waste delisting requirements may require the reporting of non-target analytes. Only after visual comparison of sample spectra with the nearest library searches may the analyst assign a tentative identification.

#### 20.6 GC/MS Tuning

Prior to any GCMS analytical sequence, including calibration, the instrument parameters for the tune and subsequent sample analyses within that sequence must be set.

Prior to tuning/auto-tuning the mass spec, the parameters may be adjusted within the specifications set by the manufacturer or the analytical method. These generally don't need any adjustment but it may be required based on the current instrument performance. If the tune verification does not pass it may be necessary to clean the source or perform additional maintenance. Any maintenance is documented in the maintenance log.

Table 20-1. Example: Instrumentation List

GC	GC/MS	Ion Trap GC/MS/MS	ICP	ICPMS	CVAA	AutoAnalyzer	IC	TOC
13	15	2	2	2	2	2	2	2

Table 20-2. Example: Software List

Lab Section	Software Title and Version
All (Data Processing and Reporting)	LIMS - TALS (TestAmerica) Ver 1.0.1.330 or higher
ECD, Fuels 8270 and VOA (Data Processing)	CHROM (TestAmerica) Ver 1.2 or higher
Data Acquisition	
ECD	Environmental ChemStation G1701DA Version D.00.01.27 Environmental ChemStation MSD ChemStation D.02.00.275
FUELS	Environmental ChemStation G1701DA Version D.00.01.27 EnviroQuant ChemStation G1701BA Version B.01.00.
8270	Environmental ChemStation MSD ChemStation D.01.00 Build 75 26-Aug-2003. Environmental ChemStation G1701DA Version D.00.01.27 08-Nov-2002 Environmental ChemStation MSD ChemStation E.01.00.237 EnviroQuant ChemStation G1702BA Version B.01.00. Varian Saturn Software Data System Finnigan Magnum Finnigan GCQ
VOA	Environmental ChemStation MSD ChemStation D.01.00.275. EnviroQuant ChemStation G1701BA Version B.01.00. Environmental ChemStation MSD ChemStation E.02.00.493. Environmental ChemStation MSD ChemStation D.02.00.275. Environmental ChemStation MSD ChemStation D.03.00.552. Environmental ChemStation G1701DA Version D.00.01.27. EnviroQuant ChemStation G1701AA Version A.03.00. Environmental ChemStation G1701CA Version C.00.00 21-Dec-1999
Metals	ELAN Version 3.0 Hotfix 3 (Build 3, 0, 6, 48d) Agilent G1834B ICP-MS ChemStation B.03.03 Perkin Elmer WinLab32 for ICP Version 4.0.0.0305 Thermo iTeVA – Analyst Version 2.2.0.51 Perkin Elmer WinLab32 for ICP Version 4.0.0.0305 WinHg Runner 1.4 CT Rev 0.286
Wet Chem	
Lachat QuikChem	Omnion Program version 3.0.220.04
Astoria Pacific	EnviroQuant ChemStation G1701BA Version B.01.00
Dionex DX-500	Dionex Peaknet Run PeakNet 5.11
Dionex ICS 2000	Chromleon Client Version 6.60 Build1428
Mitsubishi TOX-100	TOX-100 Versions 2.17 Mitsubishi Chemical Corp
OI Analytical	OI Analytical WinTOC for Solids 3.0
Teledyne Tekmar	TOCTalk Ver 4.2.105.
LECO C632	LECO C632 Carbon Determination Ver 1.60

Table 20-3. Example: Schedule of Routine Maintenance

Instrument	Procedure	Frequency
Leeman Mercury Analyzer	Check tubing for wear Fill rinse tank with 10% HCl Change drying tube Fill reductant bottle with 10% Stannous Chloride	Daily Daily As required Daily
ICP	Check pump tubing Check liquid argon supply Check fluid level in waste container Check filters Clean or replace filters Check torch Check sample spray chamber for debris Clean and align nebulizer Check entrance slit for debris Change printer ribbon Replace pump tubing	Daily Daily Daily Weekly As required Daily Monthly Monthly Monthly As required As required
ICP MS	Change pump tubing Clean torch Check / clean nebulizer Clean cones Check air filters Check multiplier voltages & do cross calibration Replace sample uptake tubing Check rotary pump oil Check oil mist filters Check chiller water level	Weekly Weekly Weekly Daily Weekly Weekly Monthly Monthly Monthly Monthly
Ion chromatograph	Check seals for leakage Replace seals/valves/lamps Replace suppressor Replace column Clean source/analyzer	Daily As required As required As required As required
UV-Vis Spectrophotometer	Clean ambient flow cell Precision check/alignment of flow cell Wavelength verification check	As required As required Semi-annually
TOC Analyzer	Check gas flow Check fluid level (IC reservoirs) Replace "O" rings Check needle Replace scrubbers (halogen and CO <sub>2</sub> ) Replace catalyst	Daily Daily As required Daily Yearly As required
Auto Analyzers	Clean sampler Check all tubing Clean inside of colorimeter Clean pump well and pump rollers Clean wash fluid receptacle Oil rollers/chains/side rails Clean optics and cells	Daily Daily Daily Quarterly Weekly Weekly Quarterly

Instrument	Procedure	Frequency
Hewlett Packard GC/MS	Ion gauge tube degassing Pump oil-level check Pump oil changing Analyzer bake-out Analyzer cleaning Resolution adjustment COMPUTER SYSTEM AND PRINTER: Air filter cleaning Change data system air filter Printer head carriage lubrication Paper sprocket cleaning Drive belt lubrication	As required Monthly Annually As required As required As required As required As required As required As required
Gas Chromatograph	Compare standard response to previous day or since last initial calibration Check carrier gas flow rate in column  Check temp. of detector, inlet, column oven Septum replacement Check system for gas leaks with SNOOP  Check for loose/frayed wires and insulation ½ Bake injector/column Change/remove sections of guard column Replace connectors/liners Change/replace column(s)	Daily  Daily via use of known compound retention Daily As required W/cylinder change as required Monthly As Required As Required As Required As Required
Electron Capture Detector (ECD)	Detector wipe test (Ni-63) Detector cleaning	Semi-annually As required
Flame Ionization Detector (FID)	Detector cleaning	As required
Balances	Class "S" traceable weight check Clean pan and check if level Field service	Daily, when used Daily At least Annually
Conductivity Meter	0.01 M KCl calibration Conductivity cell cleaning	Daily As required
Turbidimeter	Check light bulb	Daily, when used
Deionized/Distilled Water	Conductivity Point Sources Daily conductivity check Check deionizer light Monitor for VOA's System cleaning Replace cartridge & large mixed bed resins	Water Quality SOP UQA-035 Daily Daily As required As required
Drying Ovens	Temperature monitoring Temperature adjustments	Daily As required
Refrigerators/Freezers	Temperature monitoring Temperature adjustment Defrosting/cleaning	Daily As required As required
Vacuum Pumps/Air Compressor	Drained Belts checked Lubricated	Weekly Monthly Semi-annually
pH/Specific Ion Meter	Calibration/check slope Clean electrode	Daily As required

Instrument	Procedure	Frequency
BOD Incubator	Temperature monitoring Coil and incubator cleaning	Daily Monthly
Centrifuge	Check brushes and bearings	Every 6 months or as needed
Water baths	Temperature monitoring Water replaced	Daily Monthly or as needed
Zero Headspace Extractors	Verify rotation speed Check for leakage Vendor repair	As required Annually As required
TCLP Extractors	Verify rotation speed	Quarterly

## SECTION 21. MEASUREMENT TRACEABILITY

### 21.1 Overview

Traceability of measurements shall be assured using a system of documentation, calibration, and analysis of reference standards. Laboratory equipment that are peripheral to analysis and whose calibration is not necessarily documented in a test method analysis or by analysis of a reference standard shall be subject to ongoing certifications of accuracy. At a minimum, these must include procedures for checking specifications of ancillary equipment: balances, thermometers, temperature, Deionized (DI) and Reverse Osmosis (RO) water systems, automatic pipettes and other volumetric measuring devices. (Refer to Section 20.3). With the exception of Class A Glassware, quarterly accuracy checks are performed for all mechanical volumetric devices. Wherever possible, subsidiary or peripheral equipment is checked against standard equipment or standards that are traceable to national or international standards. Class A Glassware and glass microliter syringes should be routinely inspected for chips, acid etching or deformity (e.g., bent needle). If the Class A glassware or syringe is suspect, the accuracy of the glassware will be assessed prior to use.

### 21.2 NIST-Traceable Weights and Thermometers

Reference standards of measurement shall be used for calibration only and for no other purpose, unless it can be shown that their performance as reference standards would not be invalidated.

For NIST-traceable weights and thermometers, the laboratory requires that all calibrations be conducted by a calibration laboratory accredited by A2LA, NVLAP (National Voluntary Laboratory Accreditation Program), or another accreditation organization that is a signatory to a MRA (Mutual Recognition Arrangement) of one or more of the following cooperations – ILAC (International Laboratory Accreditation Cooperation) or APLAC (Asia –Pacific Laboratory Accreditation Cooperation). A calibration certificate and scope of accreditation is kept on file at the laboratory.

### 21.3 Reference Standards / Materials

Reference standards/materials, where commercially available, are traceable to certified reference materials. *Commercially prepared reference standards, to the extent available, are purchased from vendors that are accredited to ISO Guide 34 and ISO/IEC 17025. All reference standards from commercial vendors shall be accompanied with a certificate that includes at least the following information:*

- *Manufacturer*
- *Analytes or parameters calibrated*
- *Identification or lot number*
- *Calibration method*
- *Concentration with associated uncertainties*
- *Purity*

If a standard cannot be purchased from a vendor that supplies a Certificate of Analysis, the purity of the standard is documented by analysis. The receipt of all reference standards must be documented. Reference standards are labeled with a unique Standard Identification Number

## Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 7  
Effective Date: 4/23/2014  
Page 105 of 165

and expiration date. All documentation received with the reference standard is retained as a QC record and references the Standard Identification Number.

All reference, primary and working standards/materials, whether commercially purchased or laboratory prepared, must be checked regularly to ensure that the variability of the standard or material from the 'true' value does not exceed method requirements. The accuracy of calibration standards is checked by comparison with a standard from a second source. In cases where a second standard manufacturer is not available, a vendor certified different lot is acceptable for use as a second source. The appropriate Quality Control (QC) criteria for specific standards are defined in laboratory SOPs. In most cases, the analysis of an Initial Calibration Verification (ICV) or LCS (where there is no sample preparation) is used as the second source confirmation. These checks are generally performed as an integral part of the analysis method (e.g. calibration checks, laboratory control samples).

All standards and materials must be stored and handled according to method or manufacturer's requirements in order to prevent contamination or deterioration. Refer to the Corporate Environmental Health & Safety Manual or laboratory SOPs. For safety requirements, please refer to method SOPs and the laboratory Environmental Health and Safety Manual.

Standards and reference materials shall not be used after their expiration dates unless their reliability is verified by the laboratory and their use is approved by the Quality Assurance Manager. The laboratory has documented contingency procedures for re-verifying expired standards.

### 21.4 Documentation and Labeling of Standards, Reagents, and Reference Materials

Reagents must be at a minimum the purity required in the test method. The date of reagent receipt and the expiration date are documented. The lots for most of the common solvents and acids are tested for acceptability prior to company wide purchase. [Refer to TestAmerica's Corporate SOP (CA-Q-S-001), Solvent and Acid Lot Testing and Approval.]

All manufacturer or vendor supplied Certificate of Analysis or Purity must be retained, stored appropriately, and readily available for use and inspection. These records are maintained in the LIMS. Records must be kept of the date of receipt and date of expiration of standards, reagents and reference materials. In addition, records of preparation of laboratory standards, reagents, and reference materials must be retained, stored appropriately, and be readily available for use and inspection. These records are maintained in the LIMS. For detailed information on documentation and labeling, please refer to SOP TA-QA-0619, Preparation, Storage, and Verification of Standards.

Commercial materials purchased for preparation of calibration solutions, spike solutions, etc., are usually accompanied with an assay certificate or the purity is noted on the label. If the assay purity is 96% or better, the weight provided by the vendor may be used without correction. If the assay purity is less than 96% a correction will be made to concentrations applied to solutions prepared from the stock commercial material.

**21.4.1** All standards, reagents, and reference materials must be labeled in an unambiguous manner. Standards are logged into the laboratory's LIMS system, and are assigned a unique

## Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 7  
Effective Date: 4/23/2014  
Page 106 of 165

identification number. The following information is typically recorded in the electronic database within the LIMS.

- Standard ID
- Description of Standard
- Department
- Preparer's name
- Final volume and number of vials prepared
- Solvent type and lot number
- Preparation Date
- Expiration Date
- Standard source type (stock or daughter)
- Standard type (spike, surrogate, other)
- Parent standard ID (if applicable)
- Parent Standard Analyte Concentration (if applicable)
- Parent Standard Amount used (if applicable)
- Component Analytes
- Final concentration of each analyte
- Comment box (text field)

Records are maintained electronically for standard and reference material preparation. These records show the traceability to purchased stocks or neat compounds. These records also include method of preparation, date of preparation, expiration date and preparer's name or initials. Preparation procedures are provided in the Method SOPs.

**21.4.2** All standards, reagents, and reference materials must be clearly labeled with a minimum of the following information:

- Expiration Date (include prep date for reagents)
- Standard ID (from LIMS)
- Special Health/Safety warnings if applicable

Records must also be maintained of the date of receipt for commercially purchased items or date of preparation for laboratory prepared items. Special Health/Safety warnings must also be available to the analyst. This information is maintained in the LIMS and MSD Sheets.

**21.4.3** In addition, the following information may be helpful:

- Date opened (for multi-use containers, if applicable)
- Description of standard (if different from manufacturer's label or if standard was prepared in the laboratory)
- Recommended Storage Conditions
- Concentration (if applicable)

- Initials of analyst preparing standard or opening container

All containers of prepared reagents must include an expiration date and an ID number to trace back to preparation.

Procedures for preparation of reagents can be found in the Method SOPs.

Standard ID numbers must be traceable through associated logbooks, worksheets and raw data.

All reagents and standards must be stored in accordance to the following priority: 1) with the manufacturer's recommendations; 2) with requirements in the specific analytical methods as specified in laboratory SOP TA-QA-0619, Preparation, Storage, and Verification of Standards.

## SECTION 22. SAMPLING

### 22.1 Overview

The laboratory does not provide sampling services. The laboratory's responsibility in the sample collection process lies in supplying the sampler with the necessary coolers, reagent water, sample containers, preservatives, sample labels, custody seals, COC forms, ice, and packing materials required to properly preserve, pack, and ship samples to the laboratory

### 22.2 Sampling Containers

The laboratory offers clean sampling containers for use by clients. These containers are obtained from reputable container manufacturers and meet EPA specifications as required. Any certificates of cleanliness that are provided by the supplier are maintained at the laboratory.

#### 22.2.1 Preservatives

Upon request, preservatives are provided to the client in pre-cleaned sampling containers. In some cases containers may be purchased pre-preserved from the container supplier. Whether prepared by the laboratory or bought pre-preserved, the grades of the preservatives are at a minimum:

- Hydrochloric Acid – Reagent ACS (Certified VOA Free) or equivalent
- Methanol – Purge and Trap grade
- Nitric Acid – Instra-Analyzed or equivalent
- Sodium Bisulfate – ACS Grade or equivalent
- Sodium Hydroxide – Instra-Analyzed or equivalent
- Sulfuric Acid – Instra-Analyzed or equivalent
- Sodium Thiosulfate – ACS Grade or equivalent

### 22.3 Definition of Holding Time

The date and time of sampling documented on the COC form establishes the day and time zero. As a general rule, when the maximum allowable holding time is expressed in "days" (e.g., 14 days, 28 days), the holding time is based on calendar day measured. Holding times expressed in "hours" (e.g., 6 hours, 24 hours, etc.) are measured from date and time zero. The first day of holding time ends twenty-four hours after sampling. Holding times for analysis include any necessary reanalysis. However, there are some programs that determine holding time compliance based on the date and specific time of analysis compared to the time of sampling regardless of how long the holding time is.

### 22.4 Sampling Containers, Preservation Requirements, Holding Times

The preservation and holding time criteria specified in the laboratory SOPs are derived from the source documents for the methods. If method required holding times or preservation requirements are not met, the reports will be qualified using a flag, footnote or case narrative. As soon as possible or "ASAP" is an EPA designation for tests for which rapid analysis is advised, but for which neither EPA nor the laboratory have a basis for a holding time.

### 22.5 Sample Aliquots / Subsampling

Taking a representative sub-sample from a container is necessary to ensure that the analytical results are representative of the sample collected in the field. The size of the sample container, the quantity of sample fitted within the container, and the homogeneity of the sample need consideration when sub-sampling for sample preparation. It is the laboratory's responsibility to take a representative subsample or aliquot of the sample provided for analysis.

Analysts should handle each sample as if it is potentially dangerous. At a minimum, safety glasses, gloves, and lab coats must be worn when preparing aliquots for analysis.

Guidelines on taking sample aliquots & subsampling are located in SOP TA-QA-0028, Subsampling of Solid Samples.

## SECTION 23. HANDLING OF SAMPLES

Sample management procedures at the laboratory ensure that sample integrity and custody are maintained and documented from sampling/receipt through disposal.

### 23.1 Chain of Custody (COC)

The COC form is the written documented history of any sample and is initiated when bottles are sent to the field, or at the time of sampling. This form is completed by the sampling personnel and accompanies the samples to the laboratory where it is received and stored under the laboratory's custody. The purpose of the COC form is to provide a legal written record of the handling of samples from the time of collection until they are received at the laboratory. It also serves as the primary written request for analyses from the client to the laboratory. The COC form acts as a purchase order for analytical services when no other contractual agreement is in effect. An example of a COC form may be found in Figure 23-1.

### 23.1.1 Field Documentation

The information the sampler needs to provide at the time of sampling on the container label is:

- Sample identification
- Date and time
- Preservative

During the sampling process, the COC form is completed and must be legible (see Figure 23-1). This form includes information such as:

- Client name, address, phone number and fax number (if available)
- Project name and/or number
- The sample identification
- Date, time and location of sampling
- Sample collectors name
- The matrix description
- The container description
- The total number of each type of container
- Preservatives used
- Analysis requested
- Requested turnaround time (TAT)
- Any special instructions
- Purchase Order number or billing information (e.g. quote number) if available
- The date and time that each person received or relinquished the sample(s), including their signed name.

When the sampling personnel deliver the samples directly to TestAmerica personnel, the samples are stored in a cooler with ice, as applicable, and remain solely in the possession of the client's field technician until the samples are delivered to the laboratory personnel. The sample collector must assure that each container is in his/her physical possession or in his/her view at all times, or stored in such a place and manner to preclude tampering. The field technician relinquishes the samples in writing on the COC form to the sample control personnel at the laboratory or to a TestAmerica courier. When sampling personnel deliver the samples through a common carrier (Fed-Ex, UPS), the CoC relinquished date/time is completed by the field personnel and samples are released to the carrier. Samples are only considered to be received by the lab when personnel at the fixed laboratory facility have physical contact with the samples.

**Note:** Independent couriers are not required to sign the COC form. The COC is usually kept in the sealed sample cooler. The receipt from the courier is stored in log-in by date; it lists all receipts each date.

### 23.2 Sample Receipt

Samples are received at the laboratory by designated sample receiving personnel and a unique laboratory project identification number is assigned. Each sample container shall be assigned a unique sample identification number that is cross-referenced to the client identification number such that traceability of test samples is unambiguous and documented. Each sample container is affixed with a durable sample identification label. Sample acceptance, receipt, tracking and storage procedures are summarized in the following sections.


#### 23.2.1 Laboratory Receipt

When samples arrive at the laboratory, sample receiving personnel inspect the coolers and samples. The integrity of each sample must be determined by comparing sample labels or tags with the COC and by visual checks of the container for possible damage. Any non-conformance, irregularity, or compromised sample receipt must be documented on an NCM or the login Checklist in LIMS and brought to the immediate attention of the client. The COC, shipping documents, documentation of any non-conformance, irregularity, or compromised sample receipt, record of client contact, and resulting instructions become part of the project record.

##### 23.2.1.1 Unique Sample Identification

All samples that are processed through the laboratory receive a unique sample identification to ensure that there can be no confusion regarding the identity of such samples at anytime. This system includes identification for all samples, subsamples and subsequent extracts and/or digests.

The laboratory assigns a unique identification (e.g., Sample ID) code to each sample container received at the laboratory. This Primary ID is made up of the following information (consisting of 4 components):

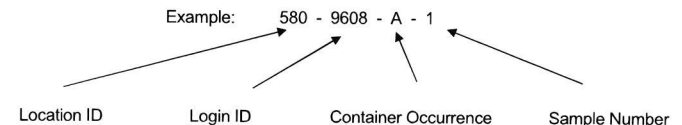

In the example above, the Location ID code identifies the lab (580 = TestAmerica Seattle), the Login ID or Job Number (unique to a particular client/job occurrence) is 9608 and the container code indicates it is the first container ("A") of Sample #1.

If the primary container goes through a prep step that creates a "new" container, then the new container is considered secondary and gets another ID. An example of this being a client sample in a 1-Liter amber bottle is sent through a Liquid/Liquid Extraction and an extraction vial is created from this step. The vial would be a SECONDARY container. The secondary ID has 5 components.

Example: 580 - 9608 - A - 1 - A ← Secondary Container Occurrence

Example: 580-9608-A-1-A, would indicate the PRIMARY container listed above that went through a step that created the 1<sup>st</sup> occurrence of a Secondary container.

With this system, a client sample can literally be tracked throughout the laboratory in every step from receipt to disposal.

### 23.3 Sample Acceptance Policy

The laboratory has a written sample acceptance policy (Figure 23-2) that clearly outlines the circumstances under which samples shall be accepted or rejected. These include:

- a COC filled out completely;
- samples must be properly labeled;
- proper sample containers with adequate volume for the analysis (Sampling Guide) and necessary QC;
- samples must be preserved according to the requirements of the requested analytical method (Sampling Guide);
- sample holding times must be adhered to (Sampling Guide);
- the project manager will be notified if any sample is received in damaged condition.

Data from samples which do not meet these criteria are flagged and the nature of the variation from policy is defined.

**23.3.1** After inspecting the samples, the sample receiving personnel sign and date the COC form, make any necessary notes of the samples' conditions and store them in appropriate refrigerators or storage locations.

**23.3.2** Any deviations from these checks that question the suitability of the sample for analysis, or incomplete documentation as to the tests required will be resolved by consultation with the client. If the sample acceptance policy criteria are not met, the laboratory shall either:

- Retain all correspondence and/or records of communications with the client regarding the disposition of rejected samples, or
- Fully document any decision to proceed with sample analysis that does not meet sample acceptance criteria.

Once sample acceptance is verified, the samples are logged into the LIMS according SOP TA-QA-0001 (Sample Receiving and Login).

### 23.4 Sample Storage

In order to avoid deterioration, contamination or damage to a sample during storage and handling, from the time of receipt until all analyses are complete, samples are stored in refrigerators, freezers or protected locations suitable for the sample matrix. Aqueous samples for metals testing are typically stored at ambient temperature. In addition, samples to be analyzed for volatile organic parameters are stored in separate refrigerators designated for

volatile organic parameters only. Samples are never to be stored with reagents, standards or materials that may create contamination.

To ensure the integrity of the samples during storage, refrigerator blanks are maintained in the volatile sample refrigerators and analyzed every two weeks. The refrigerator blanks are logged into the LIMS and treated as normal samples with the data stored and archived in the LIMS.

Analysts and technicians retrieve the sample container allocated to their analysis from the designated storage location and place them on carts, analyze the sample, and return the remaining sample or empty container to the storage location from which it originally came. All unused portions of samples, including empty sample containers, are returned to the secure sample control area. All samples are kept in the designated storage locations for two to four weeks after analysis, which meets or exceeds most sample holding times. After two to four weeks the samples are moved to dry room temperature, sample archive area where they are stored for an additional four weeks before they are disposed of. This eight week holding period allows samples to be checked if a discrepancy or question arises. Special arrangements may be made to store samples for longer periods of time. This extended holding period allows additional metal analyses to be performed on the archived sample and assists clients in dealing with legal matters or regulatory issues.

Access to the laboratory is controlled such that sample storage need not be locked at all times unless a project specifically demands it. Samples are accessible to laboratory personnel only. Visitors to the laboratory are prohibited from entering the refrigerator and laboratory areas unless accompanied by an employee of TestAmerica.

### 23.5 Hazardous Samples and Foreign Soils

To minimize exposure to personnel and to avoid potential accidents, hazardous and foreign soil samples are stored in isolated areas designated for hazardous waste or foreign soil samples only. All hazardous samples are either returned to the client or disposed of appropriately through a hazardous waste disposal firm that lab-packs all hazardous samples and removes them from the laboratory. Foreign soil samples are sent out for incineration by a USDA-approved waste disposal facility.

### 23.6 Sample Shipping

In the event that the laboratory needs to ship samples, the samples are placed in a cooler with enough ice to ensure the samples remain just above freezing and at or below 6.0°C during transit. The samples are carefully surrounded by packing material to avoid breakage (yet maintain appropriate temperature). A trip blank is enclosed for those samples requiring water/solid volatile organic analyses (see Note). The chain-of-custody form is signed by the sample control technician and attached to the shipping paperwork. Samples are generally shipped overnight express or hand-delivered by a TestAmerica courier to maintain sample integrity. All personnel involved with shipping and receiving samples must be trained to maintain the proper chain-of-custody documentation and to keep the samples intact and on ice. The Environmental, Health and Safety Manual contains additional shipping requirements.

**Note:** If a client does not request trip blank analysis on the COC or other paperwork, the laboratory will not analyze the trip blanks that were supplied. However, in the interest of good





**Figure 23-2. Example: Sample Acceptance Policy**

All incoming work will be evaluated against the criteria listed below. Where applicable, data from any samples that do not meet the criteria listed below will be noted on the laboratory report defining the nature and substance of the variation. In addition the client will be notified either by telephone, fax or e-mail ASAP after the receipt of the samples.

- 1) Samples must arrive with labels intact with a Chain of Custody filled out completely. The following information must be recorded.
  - Client name, address, phone number and fax number (if available)
  - Project name and/or number
  - The sample identification
  - Date, time and location of sampling
  - The collectors name
  - The matrix description
  - The container description
  - The total number of each type of container
  - Preservatives used
  - Analysis requested
  - Requested turnaround time (TAT)
  - Any special instructions
  - Purchase Order number or billing information (e.g. quote number) if available
  - The date and time that each person received or relinquished the sample(s), including their signed name.
  - Information must be legible
- 2) Samples must be properly labeled.
  - Use durable labels (labels provided by TestAmerica are preferred)
  - Include a unique identification number
  - Include sampling date and time & sampler ID
  - Include preservative used.
  - Use indelible ink
  - Information must be legible
- 3) Proper sample containers with adequate volume for the analysis and necessary QC are required for each analysis requested. See Lab Sampling Guide.
- 4) Samples must be preserved according to the requirements of the requested analytical method (See Sampling Guide).
- 5) Most analytical methods require chilling samples to 4° C (other than water samples for metals analysis). For these methods, the criteria are met if the samples are chilled to below 6° C and above freezing (0°C). For methods with other temperature criteria (e.g. some bacteriological methods require ≤ 10 °C), the samples must arrive within ± 2° C of the required temperature or within the method specified range. **Note:** Samples that are hand delivered to the laboratory immediately after collection may not have had time to cool sufficiently. In this case the samples will be considered acceptable as long as there is evidence that the chilling process has begun (arrival on ice).

- 5i.) Samples that are delivered to the laboratory on the same day they are collected may not meet the requirements of Section 5. In these cases, the samples shall be considered acceptable if the samples were received on ice.
  - 5ii.) If sample analysis is begun within fifteen (15) minutes of collection, thermal preservation is not required.
  - 5iii.) Thermal preservation is not required in the field if the laboratory receives and refrigerates the sample within fifteen (15) minutes of collection.
  - Chemical preservation (pH) will be verified prior to analysis and documented, either in sample control or at the analyst's level. The project manager will be notified immediately if there is a discrepancy. If analyses will still be performed, all affected results will be flagged to indicate improper preservation.
  - For Volatile Organic analyses in drinking water (Methods 502.2 or 524.2). Residual chlorine must be neutralized prior to preservation. If there is prior knowledge that the samples are not chlorinated, state it on the COC and use the VOA vials pre-preserved with HCl. The following are other options for a sampler and laboratory where the presence of chlorine is not known:
    - 1. Test for residual chlorine in the field prior to sampling.
      - If no chlorine is present, the samples are to be preserved using HCl as usual.
      - If chlorine is present, add either ascorbic acid or sodium thiosulfate prior to adding HCl.
    - 2. Use VOA vials pre-preserved with sodium thiosulfate or ascorbic acid and add HCl after filling the VOA vial with the sample.
  - **FOR WATER SAMPLES TESTED FOR CYANIDE (by Standard Methods or EPA 335)**
    - In the Field: Samples are to be tested for Sulfide using lead acetate paper prior to the addition of Sodium Hydroxide (NaOH). If sulfide is present, the sample must be treated with Cadmium Chloride and filtered prior to the addition of NaOH.
    - If the sulfide test and treatment is not performed in the field, the lab will test the samples for sulfide using lead acetate paper at the time of receipt and if sulfide is present in the sample, the client will be notified and given the option of retaking the sample and treating in the field per the method requirements or the laboratory can analyze the samples as delivered and qualify the results in the final report.
    - It is the responsibility of the client to notify the laboratory if thiosulfate, sulfite, or thiocyanate are known or suspected to be present in the sample. This notification may be on the chain of custody. The samples may need to be subcontracted to a laboratory that performs a UV digestion. If the lab does not perform the UV digestion on samples that contain these compounds, the results must be qualified in the final report.
    - The laboratory must test the sample for oxidizing agents (e.g. Chlorine) prior to analysis and treat according to the methods prior to distillation. (ascorbic acid or sodium arsenite are the preferred choice).
- 6) Sample Holding Times
    - TestAmerica will make every effort to analyze samples within the regulatory holding time. Samples must be received in the laboratory with enough time to perform the sample analysis. Except for short holding time samples (< 48hr HT), sample must be received with at least 48 hrs (working days) remaining on the holding time for us to ensure analysis.
    - Analyses that are designated as "field" analyses (Odor, pH, Dissolved Oxygen, Disinfectant Residual; a.k.a. Residual Chlorine, and Redox Potential) should be analyzed ASAP by the field sampler prior to delivering to the lab (within 15 minutes). However, if the analyses are to be

performed in the laboratory, TestAmerica will make every effort to analyze the samples within 24 hours from receipt of the samples in the testing laboratory. Samples for "field" analyses received after 4:00 pm on Friday or on the weekend will be analyzed no later than the next business day after receipt (Monday unless a holiday). Samples will remain refrigerated and sealed until the time of analysis. The actual times of all "field" sample analyses are noted on the "Short Hold Time Detail Report" in the final report. Samples analyzed in the laboratory will be qualified on the final report with an 'H' to indicate holding time exceedance.

- 7) All samples submitted for Volatile Organic analyses must have a Trip Blank submitted at the same time. TestAmerica will supply a blank with the bottle order.
- 8) The project manager will be notified if any sample is received in damaged condition. TestAmerica will request that a sample be resubmitted for analysis.
- 9) Recommendations for packing samples for shipment.
  - Pack samples in Ice rather than "Blue" ice packs.
  - Soil samples should be placed in plastic zip-lock bags. The containers often have dirt around the top and do not seal very well and are prone to intrusion from the water from melted ice.
  - Water samples would be best if wrapped with bubble-wrap or paper (newspaper, or paper towels work) and then placed in plastic zip-lock bags.
  - Fill extra cooler space with bubble wrap.

Figure 23-3. Example: Cooler Receipt Form

**Login Sample Receipt Checklist**

Client: Client X Job Number: 560-43636-1

Login Number: 43636 List Source: TestAmerica Seattle  
List Number: 1  
Creator: Blankinship, Tom X

Question	Answer	Comment
Radioactivity wasn't checked or is <1x background as measured by a survey meter.	True	
The cooler's custody seal, if present, is intact.	True	
Sample custody seals, if present, are intact.	True	
The cooler or samples do not appear to have been compromised or tampered with.	True	
Samples were received on ice.	True	
Cooler Temperature is acceptable.	True	
Cooler Temperature is recorded.	True	
COC is present.	True	
COC is filled out in ink and legible.	True	
COC is filled out with all pertinent information.	True	
Is the Field Sampler's name present on COC?	True	
There are no discrepancies between the containers received and the COC.	True	
Samples are received within Holding Time.	True	
Sample containers have legible labels.	True	
Containers are not broken or leaking.	True	
Sample collection datetimes are provided.	True	
Appropriate sample containers are used.	True	
Sample bottles are completely filled.	True	
Sample Preservation Verified.	True	
There is sufficient vol. for all requested analyses, incl. any requested MS/MS/MS.	True	
Containers requiring zero headspace have no headspace or bubble is <6mm (1/4").	N/A	
Multiphase samples are not present.	True	
Samples do not require spitting or compositing.	True	
Residual Chlorine Checked.	N/A	

TestAmerica Seattle

## SECTION 24. ASSURING THE QUALITY OF TEST RESULTS

### 24.1 Overview

In order to assure our clients of the validity of their data, the laboratory continuously evaluates the quality of the analytical process. The analytical process is controlled not only by instrument calibration as discussed in Section 20, but also by routine process quality control measurements (e.g. Blanks, Laboratory Control Samples (LCS), Matrix Spikes (MS), duplicates (DUP), surrogates, Internal Standards (IS)). These quality control checks are performed as required by the method or regulations to assess precision and accuracy. Quality control samples are to be treated in the exact same manner as the associated field samples being tested. In addition to the routine process quality control samples, Proficiency Testing (PT) Samples (concentrations unknown to laboratory) are analyzed to help ensure laboratory performance.

### 24.2 Controls

Sample preparation or pre-treatment is commonly required before analysis. Typical preparation steps include homogenization, grinding, solvent extraction, sonication, acid digestion, distillation, reflux, evaporation, drying and ashing. During these pre-treatment steps, samples are arranged into discreet manageable groups referred to as preparation (prep) batches. Prep batches provide a means to control variability in sample treatment. Control samples are added to each prep batch to monitor method performance and are processed through the entire analytical procedure with investigative/field samples.

### 24.3 Negative Controls

Table 24-1. Example – Negative Controls

Control Type	Details
Method Blank (MB)	are used to assess preparation and analysis for possible contamination during the preparation and processing steps. The specific frequency of use for method blanks during the analytical sequence is defined in the specific standard operating procedure for each analysis. Generally it is 1 for each batch of samples; not to exceed 20 environmental samples. The method blank is prepared from a clean matrix similar to that of the associated samples that is free from target analytes (e.g., Reagent water, Ottawa sand, glass beads, etc.) and is processed along with and under the same conditions as the associated samples. The method blank goes through all of the steps of the process (including as necessary: filtration, clean-ups, etc.). Reanalyze or qualify associated sample results when the concentration of a targeted analyte in the blank is at or above the reporting limit as established by the method or by regulation, AND is greater than 1/10 of the amount measured in the sample.
Calibration Blanks	are prepared and analyzed along with calibration standards where applicable. They are prepared using the same reagents that are used to prepare the standards. In some analyses the calibration blank may be included in the calibration curve.
Instrument Blanks	are blank reagents or reagent water that may be processed during an analytical sequence in order to assess contamination in the analytical system. In general, instrument blanks are used to differentiate between contamination caused by the analytical system and that caused by the sample handling or sample prep process. Instrument blanks may also be inserted throughout the analytical sequence to minimize the effect of carryover from samples with high analyte content.

Table 24-1. Example – Negative Controls

Control Type	Details
Trip Blank <sup>1</sup>	are required to be submitted by the client with each shipment of samples requiring aqueous and solid volatiles analyses (or as specified in the client's project plan). Additionally, trip blanks may be prepared and analyzed for volatile analysis of air samples, when required by the client. A trip blank may be purchased (certified clean) or is prepared by the laboratory by filling a clean container with pure deionized water that has been purged to remove any volatile compounds. Appropriate preservatives are also added to the container. The trip blank is sent with the bottle order and is intended to reflect the environment that the containers are subjected to throughout shipping and handling and help identify possible sources if contamination is found. The field sampler returns the trip blank in the cooler with the field samples.
Field Blanks <sup>1</sup>	are sometimes used for specific projects by the field samplers. A field blank prepared in the field by filling a clean container with pure reagent water and appropriate preservative, if any, for the specific sampling activity being undertaken. (EPA OSWER)
Equipment Blanks <sup>1</sup>	are also sometimes created in the field for specific projects. An equipment blank is a sample of analyte-free media which has been used to rinse common sampling equipment to check effectiveness of decontamination procedures. (TNI)
Holding Blanks	also referred to as refrigerator or freezer blanks, are used to monitor the sample storage units for volatile organic compounds during the storage of VOA samples in the laboratory

<sup>1</sup> When known, these field QC samples should not be selected for matrix QC as it does not provide information on the behavior of the target compounds in the field samples. Usually, the client sample ID will provide information to identify the field blanks with labels such as "FB", "EB", or "TB."

Evaluation criteria and corrective action for these controls are defined in the specific standard operating procedure for each analysis.

### 24.4 Positive Controls

Control samples (e.g., QC indicators) are analyzed with each batch of samples to evaluate data based upon (1) Method Performance (Laboratory Control Sample (LCS) or Blank Spike (BS)), which entails both the preparation and measurement steps; and (2) Matrix Effects (Matrix Spike (MS) or Sample Duplicate (MD, DUP), which evaluates field sampling accuracy, precision, representativeness, interferences, and the effect of the matrix on the method performed. Each regulatory program and each method within those programs specify the control samples that are prepared and/or analyzed with a specific batch.

Note that frequency of control samples vary with specific regulatory, methodology and project specific criteria. Complete details on method control samples are as listed in each analytical SOP.

#### 24.4.1 Method Performance Control - Laboratory Control Sample (LCS)

The LCS measures the accuracy of the method in a blank matrix and assesses method performance independent of potential field sample matrix effects in a laboratory batch.

The LCS is prepared from a clean matrix similar to that of the associated samples that is free from target analytes (for example: Reagent water, Ottawa sand, glass beads, etc.) and is processed along with and under the same conditions as the associated samples. The LCS is spiked with verified known amounts of analytes or is made of a material containing known and verified amounts of analytes, taken through all preparation and analysis steps along with the field samples. Where there is no preparation taken for an analysis (such as in aqueous volatiles), or when all samples and standards undergo the same preparation and analysis

process (such as Phosphorus), a calibration verification standard is reported as the LCS. In some instances where there is no practical clean solid matrix available, aqueous LCS's may be processed for solid matrices; final results may be calculated as mg/kg or ug/kg, assuming 100% solids and a weight equivalent to the aliquot used for the corresponding field samples, to facilitate comparison with the field samples.

Certified pre-made reference material purchased from a NIST/A2LA accredited vendor may also be used for the LCS when the material represents the sample matrix or the analyte is not easily spiked (e.g. solid matrix LCS for metals, TDS, etc.).

The specific frequency of use for LCS during the analytical sequence is defined in the specific standard operating procedure for each analysis. It is generally 1 for each batch of samples; not to exceed 20 environmental samples.

If the mandated or requested test method, or project requirements, do not specify the spiking components, the laboratory shall spike all reportable components to be reported in the Laboratory Control Sample (and Matrix Spike) where applicable (e.g. no spike of pH). However, in cases where the components interfere with accurate assessment (such as simultaneously spiking chlordane, toxaphene and PCBs in Method 608), the test method has an extremely long list of components or components are incompatible, at a minimum, a representative number of the listed components (see below) shall be used to control the test method. The selected components of each spiking mix shall represent all chemistries, elution patterns and masses, permit specified analytes and other client requested components. However, the laboratory shall ensure that all reported components are used in the spike mixture within a two-year time period.

- For methods that have 1-10 target analytes, spike all components.
- For methods that include 11-20 target analytes, spike at least 10 or 80%, whichever is greater.
- For methods with more than 20 target analytes, spike at least 16 components.
- Exception: Due to analyte incompatibility in pesticides, Toxaphene and Chlordane are only spiked at client request based on specific project needs.
- Exception: Due to analyte incompatibility between the various PCB Aroclors, Aroclors 1016 and 1260 are used for spiking as they cover the range of all of the Aroclors. Specific Aroclors may be used by request on a project specific basis.

## 24.5 Sample Matrix Controls

Table 24-3. Sample Matrix Control

Control Type	Details	
Matrix Spikes (MS)	Use	used to assess the effect sample matrix of the spiked sample has on the precision and accuracy of the results generated by the method used;
	Typical Frequency <sup>1</sup>	At a minimum, with each matrix-specific batch of samples processed, an MS is carried through the complete analytical procedure. Unless specified by the client, samples used for spiking are randomly selected and rotated between different client projects. If the mandated or requested test method does not specify the spiking components, the laboratory shall spike all reportable components to be reported in the Laboratory Control Sample and Matrix Spike. Refer to the method SOP for complete details
	Description	essentially a sample fortified with a known amount of the test analyte(s).
Surrogate	Use	Measures method performance to sample matrix (organics only).
	Typical Frequency <sup>1</sup>	Are added to all samples, standards, and blanks, for all organic chromatography methods except when the matrix precludes its use or when a surrogate is not available. The recovery of the surrogates is compared to the acceptance limits for the specific method. Poor surrogate recovery may indicate a problem with sample composition and shall be reported, with data qualifiers, to the client whose sample produced poor recovery.
	Description	Are similar to matrix spikes except the analytes are compounds with properties that mimic the analyte of interest and are unlikely to be found in environment samples.
Duplicates <sup>2</sup>	Use	For a measure of analytical precision, with each matrix-specific batch of samples processed, a matrix duplicate (MD or DUP) sample, matrix spike duplicate (MSD), or LCS duplicate (LCSD) is carried through the complete analytical procedure.
	Typical Frequency <sup>1</sup>	Duplicate samples are usually analyzed with methods that do not require matrix spike analysis.
	Description	Performed by analyzing two aliquots of the same field sample independently or an additional LCS.
Internal Standards	Use	Are spiked into all environmental and quality control samples (including the initial calibration standards) to monitor the qualitative aspect of organic and some inorganic analytical measurements.
	Typical Frequency <sup>1</sup>	All organic and ICP methods as required by the analytical method.
	Description	Used to correct for matrix effects and to help troubleshoot variability in analytical response and are assessed after data acquisition. Possible sources of poor internal standard response are sample matrix, poor analytical technique or instrument performance.

<sup>1</sup> See the specific analytical SOP for type and frequency of sample matrix control samples.

<sup>2</sup> LCSD's are normally not performed except when regulatory agencies or client specifications require them. The recoveries for the spiked duplicate samples must meet the same laboratory established recovery limits as the accuracy QC samples. If an LCSD is analyzed both the LCS and LCSD must meet the same recovery criteria and be included in the final report. The precision measurement is reported as "Relative Percent Difference" (RPD). Poor precision between duplicates (except LCS/LCSD) may indicate non-homogeneous matrix or sampling.

## 24.6 Acceptance Criteria (Control Limits)

As mandated by the test method and regulation, each individual analyte in the LCS, MS, or Surrogate Spike is evaluated against the control limits published in the test method. Where there are no established acceptance criteria, the laboratory calculates in-house control limits with the use of control charts or, in some cases, utilizes client project specific control limits. When this occurs, the regulatory or project limits will supersede the laboratory's in-house limits.

**Note:** For methods, analytes and matrices with very limited data (e.g., unusual matrices not analyzed often), interim limits are established using available data or by analogy to similar methods or matrices.

Once control limits have been established, they are verified, reviewed, and updated if necessary on an annual basis unless the method requires more frequent updating. Control limits are established per method (as opposed to per instrument) regardless of the number of instruments utilized.

Laboratory generated % Recovery acceptance (control) limits are generally established by taking  $\pm 3$  Standard Deviations (99% confidence level) from the average recovery of a minimum of 20-30 data points (more points are preferred).

- Regardless of the calculated limit, the limit should be no tighter than the Calibration Verification (ICV/CCV). (Unless the analytical method specifies a tighter limit).
- In-house limits cannot be any wider than those mandated in a regulated analytical method. Client or contract required control limits are evaluated against the laboratory's statistically derived control limits to determine if the data quality objectives (DQOs) can be achieved. If laboratory control limits are not consistent with DQOs, then alternatives must be considered, such as method improvements or use of an alternate analytical method.
- The lowest acceptable recovery limit will be 10% (the analyte must be detectable and identifiable). Exception: The lowest acceptable recovery limit for Benzidine will be 5% and the analyte must be detectable and identifiable.
- The maximum acceptable recovery limit will be 185%.
- The maximum acceptable RPD limit will be 35% for waters and 40% for soils. The minimum RPD limit is 10%.
- If either the high or low end of the control limit changes by  $\leq 5\%$  from previous, the control chart is visually inspected and, using professional judgment, they may be left unchanged if there is no effect on laboratory ability to meet the existing limits.

**24.6.1** The lab must be able to generate a current listing of their control limits and track when the updates are performed. In addition, the laboratory must be able to recreate historical control limits. See SOP TA-QA-0600 Quality Control Charting and Establishing Warning and Action Limits.

One example: The QA department generates a Quality Control Limit Summary that contains tables that summarize the precision and accuracy acceptability limits for analyses performed at TestAmerica Seattle. This summary includes an effective date, is updated each time new limits are generated and is located in the LIMS. Unless otherwise noted, limits within these tables are laboratory generated. The analysts are instructed to use the current limits in the laboratory (dated and approved by the Technical Director and QA Manager) and entered into the Laboratory Information Management System (LIMS). The Quality Assurance department maintains an archive of all limits used within the laboratory.

**24.6.2** A LCS that is within the acceptance criteria establishes that the analytical system is in control and is used to validate the process. Samples that are analyzed with an LCS with recoveries outside of the acceptance limits may be determined as out of control and should be reanalyzed if possible. If reanalysis is not possible, then the results for all affected analytes for samples within the same batch must be qualified when reported. The internal corrective action process (see Section 12) is also initiated if an LCS exceeds the acceptance limits. If agreed to by the client, sample results may be qualified and reported without reanalysis if:

- The analyte results are below the reporting limit and the LCS is above the upper control limit.
- If the analytical results are above the relevant regulatory limit and the LCS is below the lower control limit.

Or, for NELAC and Department Of Defense (DOD) work, there are an allowable number of Marginal Exceedances (ME):

<11 analytes	0 marginal exceedances are allowed.
11 – 30 Analytes	1 marginal exceedance is allowed
31-50 Analytes	2 marginal exceedances are allowed
51-70 Analytes	3 marginal exceedances are allowed
71-90 Analytes	4 marginal exceedances are allowed
> 90 Analytes	5 marginal exceedances are allowed

- Marginal exceedances are recovery exceedances between 3 SD and 4 SD from the mean recovery limit (TNI).
- Marginal exceedances must be random. If the same analyte exceeds the LCS control limit repeatedly, it is an indication of a systematic problem. The source of the error must be located and corrective action taken. The laboratory has a system to monitor marginal exceedances to ensure that they are random.

Though marginal exceedances may be allowed, the data must still be qualified to indicate it is outside of the normal limits.

**24.6.3** If the MS/MSDs do not meet acceptance limits, the MS/MSD and the associated spiked sample is reported with a qualifier for those analytes that do not meet limits. If obvious preparation errors are suspected, or if requested by the client, unacceptable MS/MSDs are reprocessed and reanalyzed to prove matrix interference. A more detailed discussion of acceptance criteria and corrective action can be found in the lab's method SOPs and in Section 12.

**24.6.4** If a surrogate standard falls outside the acceptance limits, if there is not obvious chromatographic matrix interference, reanalyze the sample to confirm a possible matrix effect. If the recoveries confirm or there was obvious chromatographic interference, results are reported from the original analysis and a qualifier is added. If the reanalysis meets surrogate recovery criteria, the second run is reported (or both are reported if requested by the client). Under certain circumstances, where all of the samples are from the same location and share similar chromatography, the reanalysis may be performed on a single sample rather than all of



the samples and if the surrogate meets the recovery criteria in the reanalysis, all of the affected samples would require reanalysis.

#### **24.7     Additional Procedures to Assure Quality Control**

The laboratory has written and approved method SOPs to assure the accuracy of the test method including calibration (see Section 20), use of certified reference materials (see Section 21) and use of PT samples (see Section 15).

A discussion regarding MDLs, Limit of Detection (LOD) and Limit of Quantitation (LOQ) can be found in Section 19.

- Use of formulae to reduce data is discussed in the method SOPs and in Section 20.
- Selection of appropriate reagents and standards is included in Section 9 and 21.
- A discussion on selectivity of the test is included in Section 5.
- Constant and consistent test conditions are discussed in Section 18.
- The laboratories sample acceptance policy is included in Section 23.

### **SECTION 25.     REPORTING RESULTS**

#### **25.1     Overview**

The results of each test are reported accurately, clearly, unambiguously, and objectively in accordance with State and Federal regulations as well as client requirements. Analytical results are issued in a format that is intended to satisfy customer and laboratory accreditation requirements as well as provide the end user with the information needed to properly evaluate the results. Where there is conflict between client requests and laboratory ethics or regulatory requirements, the laboratory's ethical and legal requirements are paramount, and the laboratory will work with the client during project set up to develop an acceptable solution. Refer to Section 7.

A variety of report formats are available to meet specific needs.

In cases where a client asks for simplified reports, there must be a written request from the client. There still must be enough information that would show any analyses that were out of conformance (QC out of limits) and there should be a reference to a full report that is made available to the client. Review of reported data is included in Section 19.

#### **25.2     Test Reports**

Analytical results are reported in a format that is satisfactory to the client and meets all requirements of applicable accrediting authorities and agencies. A variety of report formats are available to meet specific needs. The report is printed on laboratory letterhead, reviewed, and signed by the appropriate project manager. At a minimum, the standard laboratory report shall contain the following information:

**25.2.1**     A report title (e.g. Analytical Report For Samples) with a "sample results" column header.

**25.2.2**     Each report cover page printed on company letterhead, which includes the laboratory name, address and telephone number.

**25.2.3**     A unique identification of the report (e.g. job number) and on each page an identification in order to ensure the page is recognized as part of the report and a clear identification of the end.

**Note:**     Page numbers of report are represented as page # of ##. Where the first number is the page number and the second is the total number of pages.

**25.2.4**     A copy of the chain of custody (COC).

- Any COCs involved with Subcontracting are included.
- In most cases, the applicable COC is paginated as an integral part of the report.
- Any additional addenda to the report are also paginated as an integral part of the report (eg. Sample Receiving Checklist).

**25.2.5**     The name and address of client and a project name/number, if applicable.

**25.2.6**     Client project manager or other contact

**25.2.7**     Description and unambiguous identification of the tested sample(s) including the client identification code.

**25.2.8**     Date of receipt of sample, date and time of collection, and date(s) of test preparation and performance, and time of preparation or analysis if the required holding time for either activity is less than or equal to 72 hours.

**25.2.9**     Date reported or date of revision, if applicable.

**25.2.10**    Method of analysis including method code (EPA, Standard Methods, etc).

**25.2.11**    Reporting limit.

**25.2.12**    Method detection limits (if requested)

**25.2.13**    Definition of Data qualifiers and reporting acronyms (e.g. ND).

**25.2.14**    Sample results.

**25.2.15**    QC data consisting of method blank, surrogate, LCS, and MS/MSD recoveries and control limits.

**25.2.16** Condition of samples at receipt including temperature. This may be accomplished in a narrative or by attaching sample login sheets (Refer to Sec. 25.2.4 – Item 3 regarding additional addenda).

**25.2.17** A statement to the effect that the results relate only to the items tested and the sample as received by the laboratory.

**25.2.18** A statement that the report shall not be reproduced except in full, without prior express written approval by the laboratory coordinator.

**25.2.19** A signature and title of the person(s) accepting responsibility for the content of the report and date of issue. Signatories are appointed by the Lab Director.

**25.2.20** When NELAC accreditation is required, the lab shall certify that the test results meet all requirements of TNI or provide reasons and/or justification if they do not.

**25.2.21** Where applicable, a narrative to the report that explains the issue(s) and corrective action(s) taken in the event that a specific accreditation or certification requirement was not met.

**25.2.22** When soil samples are analyzed, a specific identification as to whether soils are reported on a "wet weight" or "dry weight" basis.

**25.2.23** Appropriate laboratory certification number for the state of origin of the sample, if applicable.

**25.2.24** If only part of the report is provided to the client (client requests some results before all of it is complete), it must be clearly indicated on the report (e.g., preliminary). A complete report must be sent once all of the work has been completed.

**25.2.25** Any non-TestAmerica subcontracted analysis results are provided as a separate report on the official letterhead of the subcontractor. All TestAmerica subcontracting is clearly identified on the report as to which laboratory performed a specific analysis.

**25.2.26** A clear statement notifying the client that non-accredited tests were performed and directing the client to the laboratory's accreditation certificates of approval shall be provided when non-accredited tests are included in the report.

Note: Refer to the Corporate SOP on Electronic Reporting and Signature Policy (No. CA-I-P-002) for details on internally applying electronic signatures of approval.

## **25.3 Reporting Level or Report Type**

The laboratory offers four levels of quality control reporting. Each level, in addition to its own specific requirements, contains all the information provided in the preceding level. The packages provide the following information in addition to the information described above:

- Level I is a report with the features described in Section 25.2 above.
- Level II is a Level I report plus summary information, including results for the method blank reported to the laboratory MDL, percent recovery for laboratory control samples and matrix spike samples, and the RPD values for all MSD and sample duplicate analyses.
- Level III contains all the information supplied in Level II, but presented on the CLP-like summary forms, and relevant calibration information. A Level II report is not included, unless specifically requested. No raw data is provided.
- Level IV is the same as Level III with the addition of all raw supporting data.

In addition to the various levels of QC packaging, the laboratory also provides reports in diskette deliverable form. Initial reports may be provided to clients by facsimile. All faxed reports are followed by hardcopy. Procedures used to ensure client confidentiality are outlined in Section 25.6.

### **25.3.1 Electronic Data Deliverables (EDDs)**

EDDs are routinely offered as part of TestAmerica's services. TestAmerica Seattle offers a variety of EDD formats including Environmental Restoration Information Management System (ERPIMS), Staged Electronic Data Deliverable (SEDD) Environmental Quality Information System (EQulS), Electronic Deliverable Format (EDF), Excel and custom files.

EDD specifications are submitted to the IT department by the PM for review and undergo the contract review process. Once the facility has committed to providing data in a specific electronic format, the coding of the format may need to be performed. This coding is documented and validated. The validation of the code is retained by the IT staff coding the EDD.

EDDs shall be subject to a review to ensure their accuracy and completeness. If EDD generation is automated, review may be reduced to periodic screening if the laboratory can demonstrate that it can routinely generate that EDD without errors. Any revisions to the EDD format must be reviewed until it is demonstrated that it can routinely be generated without errors. If the EDD can be reproduced accurately and if all subsequent EDDs can be produced error-free, each EDD does not necessarily require a review.

## **25.4 Supplemental Information for Test**

The lab identifies any unacceptable QC analyses or any other unusual circumstances or observations such as environmental conditions and any non-standard conditions that may have affected the quality of a result. This is typically in the form of a footnote or a qualifier and/or a narrative explaining the discrepancy in the front of the report.

Numeric results with values outside of the calibration range, either high or low are qualified as 'estimated'.

Where quality system requirements are not met, a statement of compliance/non-compliance with requirements and/or specifications is required, including identification of test results derived from any sample that did not meet TNI sample acceptance requirements such as improper container, holding time, or temperature.

Where applicable, a statement on the estimated uncertainty of measurements; information on uncertainty is needed when a client's instructions so require.

Opinions and Interpretations - The test report contains objective information, and generally does not contain subjective information such as opinions and interpretations. If such information is required by the client, the Laboratory Director will determine if a response can be prepared. If so, the Laboratory Director will designate the appropriate member of the management team to prepare a response. The response will be fully documented, and reviewed by the Laboratory Director, before release to the client. There may be additional fees charged to the client at this time, as this is a non-routine function of the laboratory.

**Note:** Review of data deliverable packages for submittal to regulatory authorities requires responses to non-conforming data concerning potential impact on data quality. This necessitates a limited scope of interpretation, and this work is performed by the QA Department. This is the only form of "interpretation" of data that is routinely performed by the laboratory.

When opinions or interpretations are included in the report, the laboratory provides an explanation as to the basis upon which the opinions and interpretations have been made. Opinions and interpretations are clearly noted as such and where applicable, a comment should be added suggesting that the client verify the opinion or interpretation with their regulator.

#### **25.5 Environmental Testing Obtained From Subcontractors**

If the laboratory is not able to provide the client the requested analysis, the samples would be subcontracted following the procedures outlined in the Corporate SOP on Subcontracting (SOP No. CA-L-S-002).

Data reported from analyses performed by a subcontractor laboratory are clearly identified as such on the analytical report provided to the client. Results from a subcontract laboratory outside of TestAmerica are reported to the client on the subcontract laboratory's original report stationary and the report includes any accompanying documentation.

#### **25.6 Client Confidentiality**

In situations involving the transmission of environmental test results by telephone, facsimile or other electronic means, client confidentiality must be maintained.

TestAmerica will not intentionally divulge to any person (other than the Client or any other person designated by the Client in writing) any information regarding the services provided by TestAmerica or any information disclosed to TestAmerica by the Client. Furthermore,

information known to be potentially endangering to national security or an entity's proprietary rights will not be released.

**Note:** This shall not apply to the extent that the information is required to be disclosed by TestAmerica under the compulsion of legal process. TestAmerica will, to the extent feasible, provide reasonable notice to the client before disclosing the information.

**Note:** Authorized representatives of an accrediting authority are permitted to make copies of any analyses or records relevant to the accreditation process, and copies may be removed from the laboratory for purposes of assessment.

**25.6.1** Report deliverable formats are discussed with each new client. If a client requests that reports be faxed or e-mailed, the reports are faxed with a cover sheet or e-mailed with the following note that includes a confidentiality statement similar to the following:

This report is issued solely for the use of the person or company to whom it is addressed. Any use, copying or disclosure other than by the intended recipient is unauthorized. If you have received this report in error, please notify the sender immediately at 253-922-2310 and destroy this report immediately.

#### **25.7 Format of Reports**

The format of reports is designed to accommodate each type of environmental test carried out and to minimize the possibility of misunderstanding or misuse.

#### **25.8 Amendments to Test Reports**

Corrections, additions, or deletions to reports are only made when justification arises through supplemental documentation. Justification is documented using the laboratory's corrective action system (refer to Section 12).

The revised report is retained on the Archive data server, as is the original report. The revised report is stored in the Archive data server under the job number followed by the appropriate revision number. The revised report will have the word "revised" or "amended" next to the date rather than the word "reported".

When the report is re-issued, a notation of "Revision" with the appropriate number is placed on the cover/signature page of the report or at the top of the narrative page with a brief explanation of reason for the re-issue and a reference back to the last final report generated. For example: Report was revised on 11/3/08 to include toluene in sample NQA1504 per client's request. This final report replaces the final report generated on 10/27/08 at 10:47am.

#### **25.9 Policies on Client Requests for Amendments**

##### **25.9.1 Policy on Data Omissions or Reporting Limit Increases**


Fundamentally, our policy is simply to not omit previously reported results (including data qualifiers) or to not raise reporting limits and report sample results as ND. This policy has few exceptions. Exceptions are:

- Laboratory error.
- Sample identification is indeterminate (confusion between COC and sample labels).
- An incorrect analysis (not analyte) was requested (e.g., COC lists 8315 but client wanted 8310). A written request for the change is required.
- Incorrect limits reported based on regulatory requirements.
- The requested change has absolutely no possible impact on the interpretation of the analytical results and there is no possibility of the change being interpreted as misrepresentation by anyone inside or outside of our company.

#### 25.9.2 Multiple Reports

TestAmerica does not issue multiple reports for the same work order where there is different information on each report (this does not refer to copies of the same report) unless required to meet regulatory needs and approved by QA.

#### Appendix 1. Laboratory Floor Plan

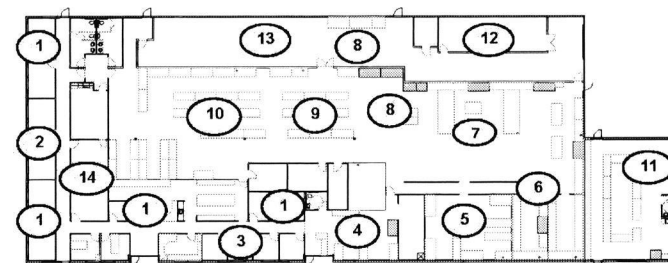

Key Areas *		
1. Admin/HR	6. Wet / General Chem	11. GC & GC-MS VOA
2. QA Offices	7. Metals Prep	12. Sample Disposal
3. Project Management	8. Organics Prep	13. GeoTech. Storage, Warehouse
4. Sample Control	9. GC-MS SVOA	14. Conference Room
5. Metals & Mercury Lab	10. GC SVOA	

**Appendix 2. Glossary/Acronyms (EL-V1M2 Sec. 3.1)****Glossary:**

**Acceptance Criteria:** Specified limits placed on characteristics of an item, process, or service defined in requirement documents. (ASQC)

**Accreditation:** The process by which an agency or organization evaluates and recognizes a laboratory as meeting certain predetermined qualifications or standards, thereby accrediting the laboratory.

**Accuracy:** The degree of agreement between an observed value and an accepted reference value. Accuracy includes a combination of random error (precision) and systematic error (bias) components which are due to sampling and analytical operations; a data quality indicator. (QAMS)

**Analyst:** The designated individual who performs the "hands-on" analytical methods and associated techniques and who is the one responsible for applying required laboratory practices and other pertinent quality controls to meet the required level of quality.

**Analytical Uncertainty:** A subset of Measurement Uncertainty that includes all laboratory activities performed as part of the analysis. (TNI)

**Assessment:** The evaluation process used to measure or establish the performance, effectiveness, and conformance of an organization and/or its systems to defined criteria (to the standards and requirements of laboratory accreditation). (TNI)

**Audit:** A systematic and independent examination of facilities, equipment, personnel, training, procedures, record-keeping, data validation, data management, and reporting aspects of a system to determine whether QA/QC and technical activities are being conducted as planned and whether these activities will effectively achieve quality objectives. (TNI)

**Batch:** Environmental samples that are prepared and/or analyzed together with the same process and personnel, using the same lot(s) of reagents. A **preparation batch** is composed of one (1) to twenty (20) environmental samples of the same quality systems matrix, meeting the above mentioned criteria and with a maximum time between the start of processing of the first and last sample in the batch to be twenty-four (24) hours. An **analytical batch** is composed of prepared environmental samples (extracts, digestates or concentrates) which are analyzed together as a group. An analytical batch can include prepared samples originating from various quality system matrices and can exceed twenty (20) samples. (TNI)

**Bias:** The systematic or persistent distortion of a measurement process, which causes errors in one direction (i.e., the expected sample measurement is different from the sample's true value). (TNI)

**Blank:** A sample that has not been exposed to the analyzed sample stream in order to monitor contamination during sampling, transport, storage or analysis. The blank is subjected to the usual analytical and measurement process to establish a zero baseline or background value and is sometimes used to adjust or correct routine analytical results. (ASQC)

**Calibration:** A set of operations that establish, under specified conditions, the relationship between values of quantities indicated by a measuring instrument or measuring system, or values represented by a material measure or a reference material, and the corresponding values realized by standards. (TNI)

1) In calibration of support equipment the values realized by standards are established through the use of reference standards that are traceable to the International System of Units (SI).

2) In calibration according to methods, the values realized by standards are typically established through the use of Reference Materials that are either purchased by the laboratory with a certificate of analysis or purity, or prepared by the laboratory using support equipment that has been calibrated or verified to meet specifications.

**Calibration Curve:** The mathematical relationship between the known values, such as concentrations, of a series of calibration standards and their instrument response. (TNI)

**Calibration Standard:** A substance or reference material used to calibrate an instrument (QAMS)

**Certified Reference Material (CRM):** A reference material accompanied by a certificate, having a value, measurement uncertainty, and stated metrological traceability chain to a national metrology institute. (TNI)

**Chain of Custody (COC) Form:** Record that documents the possession of the samples from the time of collection to receipt in the laboratory. This record generally includes: the number and types of containers; the mode of collection; the collector; time of collection; preservation; and requested analyses. (TNI)

**Compromised Samples:** Those samples which are improperly sampled, insufficiently documented (chain of custody and other sample records and/or labels), improperly preserved, collected in improper containers, or exceeding holding times when delivered to a laboratory. Under normal conditions, compromised samples are not analyzed. If emergency situation require analysis, the results must be appropriately qualified.

**Confidential Business Information (CBI):** Information that an organization designates as having the potential of providing a competitor with inappropriate insight into its management, operation or products. TNI and its representatives agree to safeguarding identified CBI and to maintain all information identified as such in full confidentiality.

**Confirmation:** Verification of the identity of a component through the use of an approach with a different scientific principle from the original method. These may include, but are not limited to Second Column Confirmation; Alternate wavelength; Derivatization; Mass spectral interpretation; Alternative detectors or Additional Cleanup procedures. (TNI)

**Conformance:** An affirmative indication or judgment that a product or service has met the requirements of the relevant specifications, contract, or regulation; also the state of meeting the requirements. (ANSI/ASQC E4-1994)

**Correction:** Actions necessary to correct or repair analysis specific non-conformances. The acceptance criteria for method specific QC and protocols as well as the associated corrective actions. The analyst will most frequently be the one to identify the need for this action as a result of calibration checks and QC sample analysis. No significant action is taken to change behavior, process or procedure.

**Corrective Action:** The action taken to eliminate the causes of an existing nonconformity, defect or other undesirable situation in order to prevent recurrence. (ISO 8402)

**Data Audit:** A qualitative and quantitative evaluation of the documentation and procedures associated with environmental measurements to verify that the resulting data re of acceptable quality (i.e., that they meet specified acceptance criteria).

**Data Reduction:** The process of transforming the number of data items by arithmetic or statistical calculations, standard curves, and concentration factors, and collation into a more useable form. (TNI)

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 7  
Effective Date: 4/23/2014  
Page 135 of 165

**Deficiency:** An unauthorized deviation from acceptable procedures or practices, or a defect in an item. (ASQC)

**Demonstration of Capability:** A procedure to establish the ability of the analyst to generate analytical results of acceptable accuracy and precision. (TNI)

**Document Control:** The act of ensuring that documents (and revisions thereto) are proposed, reviewed for accuracy, approved for release by authorized personnel, distributed properly, and controlled to ensure use of the correct version at the location where the prescribed activity is performed. (ASQC)

**Duplicate Analyses:** The analyses or measurements of the variable of interest performed identically on two subsamples of the same sample. The results from duplicate analyses are used to evaluate analytical or measurement precision but not the precision of sampling, preservation or storage internal to the laboratory. (EPA-QAD)

**Equipment Blank:** Sample of analyte-free media which has been used to rinse common sampling equipment to check effectiveness of decontamination procedures.

**External Standard Calibration:** Calibrations for methods that do not utilize internal standards to compensate for changes in instrument conditions.

**Field Blank:** Blank prepared in the field by filling a clean container with pure de-ionized water and appropriate preservative, if any, for the specific sampling activity being undertaken (EPA OSWER)

**Field of Accreditation:** Those matrix, technology/method, and analyte combinations for which the accreditation body offers accreditation.

**Holding Times:** The maximum time that samples may be held prior to analyses and still be considered valid or not compromised. (40 CFR Part 136)

**Internal Standard:** A known amount of standard added to a test portion of a sample as a reference for evaluating and controlling the precision and bias of the applied analytical test method. (TNI)

**Internal Standard Calibration:** Calibrations for methods that utilize internal standards to compensate for changes in instrument conditions.

**Instrument Blank:** A clean sample (e.g., distilled water) processed through the instrumental steps of the measurement process; used to determine instrument contamination. (EPA-QAD)

**Instrument Detection Limit (IDL):** The minimum amount of a substance that can be measured with a specified degree of confidence that the amount is greater than zero using a specific instrument. The IDL is associated with the instrumental portion of a specific method only, and sample preparation steps are not considered in its derivation. The IDL is a statistical estimation at a specified confidence interval of the concentration at which the relative uncertainty is  $\pm 100\%$ . The IDL represents a range where qualitative detection occurs on a specific instrument. Quantitative results are not produced in this range.

**Laboratory Control Sample** (however named, such as laboratory fortified blank, spiked blank, or QC check sample): A sample matrix, free from the analytes of interest, spiked with verified known amounts of analytes or a material containing known and verified amounts of analytes, taken through all preparation and analysis steps of the procedure unless otherwise noted in a reference method. It is generally used to establish intra-laboratory or analyst specific precision and bias or to assess the performance of all or a portion of the measurement system.

Company Confidential & Proprietary

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 7  
Effective Date: 4/23/2014  
Page 136 of 165

An LCS shall be prepared at a minimum of 1 per batch of 20 or less samples per matrix type per sample extraction or preparation method except for analytes for which spiking solutions are not available such as total suspended solids, total dissolved solids, total volatile solids, total solids, pH, color, odor, temperature, dissolved oxygen or turbidity. The results of these samples shall be used to determine batch acceptance.

**Least Squares Regression (1<sup>st</sup> Order Curve):** The least squares regression is a mathematical calculation of a straight line over two axes. The y axis represents the instrument response (or Response ratio) of a standard or sample and the x axis represents the concentration. The regression calculation will generate a correlation coefficient (r) that is a measure of the "goodness of fit" of the regression line to the data. A value of 1.00 indicates a perfect fit. In order to be used for quantitative purposes, r must be greater than or equal to 0.99 for organics and 0.995 for inorganics.

**Limit(s) of Detection (LOD) [a.k.a., Method Detection Limit (MDL)]:** A laboratory's estimate of the minimum amount of an analyte in a given matrix that an analytical process can reliably detect in their facility. (TNI)

**LOD Verification [a.k.a., MDL Verification]:** A processed QC sample in the matrix of interest, spiked with the analyte at no more than 3X the LOD for single analyte tests and 4X the LOD for multiple analyte tests and processed through the entire analytical procedure.

**Limit(s) of Quantitation (LOQ) [a.k.a., Reporting Limit]:** The minimum levels, concentrations, or quantities of a target variable (e.g., target analyte) that can be reported with a specified degree of confidence. (TNI)

**(QS) Matrix:** The component or substrate that contains the analyte of interest. For purposes of batch and QC requirement determinations, the following matrix distinctions shall be used:

**Aqueous:** Any aqueous sample excluded from the definition of Drinking Water or Saline/Estuarine. Includes surface water, groundwater, effluents, and TCLP or other extracts.

**Drinking Water:** Any aqueous sample that has been designated as a potable or potential potable water source.

**Saline/Estuarine:** Any aqueous sample from an ocean or estuary, or other salt water source such as the Great Salt Lake.

**Non-Aqueous Liquid:** Any organic liquid with <15% settleable solids.

**Biological Tissue:** Any sample of a biological origin such as fish tissue, shellfish, or plant material. Such samples shall be grouped according to origin.

**Solids:** Includes soils, sediments, sludges, and other matrices with >15% settleable solids.

**Chemical Waste:** A product or by-product of an industrial process that results in a matrix not previously defined.

**Air & Emissions:** Whole gas or vapor samples including those contained in flexible or rigid wall containers and the extracted concentrated analytes of interest from a gas or vapor that are collected with a sorbent tube, impinger solution, filter, or other device. (TNI)

**Matrix Spike (spiked sample or fortified sample):** A sample prepared, taken through all sample preparation and analytical steps of the procedure unless otherwise noted in a referenced method, by adding a known amount of target analyte to a specified amount of sample for which an independent test

Company Confidential & Proprietary



## Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 7  
Effective Date: 4/23/2014  
Page 137 of 165

result of target analyte concentration is available. Matrix spikes are used, for example, to determine the effect of the matrix on a method's recovery efficiency.

**Matrix Spike Duplicate (spiked sample or fortified sample duplicate):** A replicate matrix spike is prepared and analyzed to obtain a measure of the precision of the recovery for each analyte.

**Method Blank:** A sample of a matrix similar to the batch of associated samples (when available) that is free from the analytes of interest and is processed simultaneously with and under the same conditions as samples through all steps of the analytical procedures, and in which no target analytes or interferences are present at concentrations that impact the analytical results for sample analyses.

**Method Detection Limit:** The minimum concentration of a substance (an analyte) that can be measured and reported with 99% confidence that the analyte concentration is greater than zero and is determined from analysis of a sample in a given matrix containing the analyte. (40 CFR Part 136, Appendix B)

**Negative Control:** Measures taken to ensure that a test, its components, or the environment do not cause undesired effects, or produce incorrect test results.

**Non-conformance:** An indication, judgment, or state of not having met the requirements of the relevant specifications, contract, or regulation.

**Performance Audit:** The routine comparison of independently obtained qualitative and quantitative measurement system data with routinely obtained data in order to evaluate the proficiency of an analyst or laboratory.

**Positive Control:** Measures taken to ensure that a test and/or its components are working properly and producing correct or expected results from positive test subjects.

**Precision:** The degree to which a set of observations or measurements of the same property, obtained under similar conditions, conform to themselves; a data quality indicator. Precision is usually expressed as standard deviation, variance or range, in either absolute or relative terms. (TNI)

**Preservation:** Any conditions under which a sample must be kept in order to maintain chemical and/or biological integrity prior to analysis. (TNI)

**Proficiency Testing:** A means of evaluating a laboratory's performance under controlled conditions relative to a given set of criteria through analysis of unknown samples provided by an external source. (TNI)

**Proficiency Testing Program:** The aggregate of providing rigorously controlled and standardized environmental samples to a laboratory for analysis, reporting of results, statistical evaluation of the results and the collective demographics and results summary of all participating laboratories. (TNI)

**Proficiency Test Sample (PT):** A sample, the composition of which is unknown to the laboratory and is provided to test whether the laboratory can produce analytical results within specified acceptance criteria. (TNI)

**Quality Assurance:** An integrated system of management activities involving planning, implementation, assessment, reporting and quality improvement to ensure that a process, item, or service is of the type of quality needed and expected by the client. (TNI)

**Quality Assurance [Project] Plan (QAPP):** A formal document describing the detailed quality control procedures by which the quality requirements defined for the data and decisions pertaining to a specific project are to be achieved. (EAP-QAD)

Company Confidential & Proprietary

## Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 7  
Effective Date: 4/23/2014  
Page 138 of 165

**Quality Control:** The overall system of technical activities that measures the attributes and performance of a process, item, or service against defined standards to verify that they meet the stated requirements established by the customer; operational techniques and activities that are used to fulfill requirements for quality; also the system of activities and checks used to ensure that measurement systems are maintained within prescribed limits, providing protection against "out of control" conditions and ensuring that the results are of acceptable quality. (TNI)

**Quality Control Sample:** A sample used to assess the performance of all or a portion of the measurement system. One of any number of samples, such as Certified Reference Materials, a quality system matrix fortified by spiking, or actual samples fortified by spiking, intended to demonstrate that a measurement system or activity is in control. (TNI)

**Quality Manual:** A document stating the management policies, objectives, principles, organizational structure and authority, responsibilities, accountability, and implementation of an agency, organization, or laboratory, to ensure the quality of its product and the utility of its product to its users. (TNI)

**Quality System:** A structured and documented management system describing the policies, objectives, principles, organizational authority, responsibilities, accountability, and implementation plan of an organization for ensuring quality in its work processes, products (items), and services. The quality system provides the framework for planning, implementing, and assessing work performed by the organization and for carrying out required QA and QC activities. (TNI)

**Raw Data:** The documentation generated during sampling and analysis. This documentation includes, but is not limited to, field notes, electronic data, magnetic tapes, untabulated sample results, QC sample results, print outs of chromatograms, instrument outputs, and handwritten records. (TNI)

**Record Retention:** The systematic collection, indexing and storing of documented information under secure conditions.

**Reference Material:** Material or substance, one or more properties of which are sufficiently homogeneous and well established to be used for the calibration of an apparatus, the assessment of a measurement method, or for assigning values to materials. (TNI)

**Reference Standard:** Standard used for the calibration of working measurement standards in a given organization or a given location. (TNI)

**Sampling:** Activity related to obtaining a representative sample of the object of conformity assessment, according to a procedure.

**Second Order Polynomial Curve (Quadratic):** The 2<sup>nd</sup> order curves are a mathematical calculation of a slightly curved line over two axis. The y axis represents the instrument response (or Response ratio) of a standard or sample and the x axis represents the concentration. The 2<sup>nd</sup> order regression will generate a coefficient of determination (COD or  $r^2$ ) that is a measure of the "goodness of fit" of the quadratic curvature the data. A value of 1.00 indicates a perfect fit. In order to be used for quantitative purposes,  $r^2$  must be greater than or equal to 0.99.

**Selectivity:** The ability to analyze, distinguish, and determine a specific analyte or parameter from another component that may be a potential interferent or that may behave similarly to the target analyte or parameter within the measurement system. (TNI)

**Sensitivity:** The capability of a method or instrument to discriminate between measurement responses representing different levels (e.g., concentrations) of a variable of interest. (TNI)

Company Confidential & Proprietary

## Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 7  
Effective Date: 4/23/2014  
Page 139 of 165

**Spike:** A known mass of target analyte added to a blank, sample or sub-sample; used to determine recovery efficiency or for other quality control purposes.

**Standard:** The document describing the elements of laboratory accreditation that has been developed and established within the consensus principles of standard setting and meets the approval requirements of standard adoption organizations procedures and policies. (TNI)

**Standard Operating Procedures (SOPs):** A written document which details the method for an operation, analysis, or action, with thoroughly prescribed techniques and steps. SOPs are officially approved as the methods for performing certain routine or repetitive tasks. (TNI)

**Storage Blank:** A blank matrix stored with field samples of a similar matrix (volatiles only) that measures storage contribution to any source of contamination.

**Surrogate:** A substance with properties that mimic the analyte of interest. It is unlikely to be found in environment samples and is added to them for quality control purposes.

Surrogate compounds must be added to all samples, standards, and blanks, for all organic chromatography methods except when the matrix precludes its use or when a surrogate is not available. Poor surrogate recovery may indicate a problem with sample composition and shall be reported to the client whose sample produced poor recovery. (QAMS)

**Systems Audit (also Technical Systems Audit):** A thorough, systematic, qualitative on-site assessment of the facilities, equipment, personnel, training, procedures, record keeping, data validation, data management, and reporting aspects of a total measurement system. (EPA-QAD)

**Technical Manager:** A member of the staff of an environmental laboratory who exercises actual day-to-day supervision of laboratory operations for the appropriate fields of accreditation and reporting of results

**Technology:** A specific arrangement of analytical instruments, detection systems, and/or preparation techniques.

**Traceability:** The ability to trace the history, application, or location of an entity by means of recorded identifications. In a calibration sense, traceability relates measuring equipment to national or international standards, primary standards, basic physical constants or properties, or reference materials. In a data collection sense, it relates calculations and data generated throughout the project back to the requirements for the quality of the project. (TNI)

**Trip Blank:** A blank matrix placed in a sealed container at the laboratory that is shipped, held unopened in the field, and returned to the laboratory in the shipping container with the field samples.

**Uncertainty:** A parameter associated with the result of a measurement that characterizes the dispersion of the value that could reasonably be attributed to the measured value.

### Acronyms:

CAR – Corrective Action Report  
CCV – Continuing Calibration Verification  
CF – Calibration Factor  
CFR – Code of Federal Regulations  
COC – Chain of Custody  
DOC – Demonstration of Capability  
DQO – Data Quality Objectives

Company Confidential & Proprietary

## Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 7  
Effective Date: 4/23/2014  
Page 140 of 165

DUP - Duplicate  
EHS – Environment, Health and Safety  
EPA – Environmental Protection Agency  
GC - Gas Chromatography  
GC/MS - Gas Chromatography/Mass Spectrometry  
HPLC - High Performance Liquid Chromatography  
ICP - Inductively Coupled Plasma Atomic Emission Spectroscopy  
ICP/MS – ICP/Mass Spectrometry  
ICV – Initial Calibration Verification  
IDL – Instrument Detection Limit  
IH – Industrial Hygiene  
IS – Internal Standard  
LCS – Laboratory Control Sample  
LCSD – Laboratory Control Sample Duplicate  
LIMS – Laboratory Information Management System  
LOD – Limit of Detection  
LOQ – Limit of Quantitation  
LOQV - Limit of Quantitation Check Standard  
MDL – Method Detection Limit  
MDLV – MDL Verification Check Standard  
MRL – Method Reporting Limit  
MS – Matrix Spike  
MSD – Matrix Spike Duplicate  
MSDS - Material Safety Data Sheet  
NELAP - National Environmental Laboratory Accreditation Program  
PT – Performance Testing  
TNI – The NELAC Institute  
QAM – Quality Assurance Manual  
QA/QC – Quality Assurance / Quality Control  
QAPP – Quality Assurance Project Plan  
RF – Response Factor  
RPD – Relative Percent Difference  
RSD – Relative Standard Deviation  
SD – Standard Deviation  
SOP – Standard Operating Procedure  
TAT – Turn-Around-Time  
VOA – Volatiles  
VOC – Volatile Organic Compound

Company Confidential & Proprietary

### Appendix 3. Laboratory Certifications, Accreditations, Validations

TestAmerica Seattle maintains accreditations, certifications, and approvals with numerous state and national entities. Programs vary but may include on-site audits, reciprocal agreements with another entity, performance testing evaluations, review of the QA Manual, Standard Operating Procedures, Method Detection Limits, training records, etc. At the time of this QA Manual revision, the laboratory has accreditation/certification/licensing with the following organizations:

Organization	Lab ID Number
DoD ELAP	L2236
ISO 17025	L2236
Alaska	UST-022
California (NELAP)	01115CA
Montana	(UST – no number)
Oregon (NELAP)	WA100007
Washington	C553
USDA Soil Permit	P330-11-00222
USFWS Tissue Import Permit	LE192332-0

The certificates and parameter lists (which may differ) for each organization are available, upon request, from a laboratory representative or may be found on the corporate web site, the laboratory's public server, the final report review table, and in the following offices: QA, marketing, and project management.

### Appendix 4: Summary of Calibration, QC and Corrective Action Procedures for GC Organics

Method	QC Check	Frequency	Acceptance Criteria <sup>1</sup>	Corrective Action <sup>1</sup>
SW8081A SW8082	Minimum five-point initial calibration for all target analytes <sup>2</sup>	Initial calibration prior to sample analysis. Perform instrument re-calibration once per year minimum.	Linear regression correlation coefficient: $r \geq 0.990$ , $r^2 \geq 0.990$ . RSD of CF $\leq 20\%$ . <b>For DoD:</b> Linear regression correlation coefficient: $r \geq 0.995$ , $r^2 \geq 0.990$ . RSD of CF $\leq 20\%$ .	Correct problem then repeat initial calibration
	Initial calibration verification (ICV) must be from a 2 <sup>nd</sup> source	Once immediately following initial calibration	All target analytes within 15% of expected value. <b>For DoD:</b> All target analytes within 20% of expected value.	Correct problem then repeat initial calibration
	Continuing calibration verification (CCV)	Before sample analysis, after every 10 samples, and at the end of the analysis sequence	All analytes within 15% of expected value and within the RT Window. <b>For DoD:</b> All target analytes within 20% of expected value.	Correct problem then repeat initial CCV (re-calibrate if necessary) and re-analyze all samples since last successful CCV.
	Breakdown check (Endrin and DDT) <sup>3</sup>	Before sample analysis	Degradation $\leq 15\%$ for either Endrin or DDT.	Inlet/column maintenance; repeat breakdown check and re-analyze all samples since last successful breakdown check.
	Method blank	One per analytical prep batch, not to exceed 20 samples in a batch.	No analytes detected $\geq \frac{1}{2}$ RL or MDL, whichever is greater <sup>4</sup> .	Correct problem then re-prepare and analyze method blank and all samples processed with the contaminated blank.
	LCS for all analytes, must be from a 2 <sup>nd</sup> source.	One per prep batch, not to exceed 20 samples in a batch.	See Control Limits	Re-prepare and analyze the LCS and all samples in the affected analytical batch.
	Surrogate(s)	Every sample, spike, standard, and method blank.	See Control Limits	Check system, re-inject, re-extract <sup>5</sup> .
	MS/MSD, must be from a 2 <sup>nd</sup> source. Rotate Aroclors each quarter.	One per batch per matrix, if insufficient sample for MS/MSD, then a LCS/LCSD will be analyzed.	See Control Limits	None (LCS is used to determine if data is acceptable).
	Second-column confirmation	100% for all positive results	Same as for initial or primary column analysis	Same as for initial or primary column analysis. If the relative % difference of results between the 2 columns is greater than 40%, a comment should be placed in LIMS.
	Retention time window calculated for each analyte (see section 9 for how to calculate RTWs).	System set-up, with each new column or major instrument maintenance. Update the mid-RTW at the start of the run or daily.	Each analyte of the LCS, MS/MSD and CCV must be within the calculated RTW.	Correct the problem and re-process or re-analyze samples. If questions, see the supervisor or Technical Manager.
	MDL verification	Minimum yearly <b>For DOD:</b> Minimum quarterly	Detectable	Re-evaluate MDL standard used and MDL; see Technical Manager.
	LOQ verification	<b>For DOD:</b> Minimum quarterly	Detectable	Re-evaluate LOQ standard used and LOQ; see Technical Manager.

<sup>1</sup> – 8081A only

Company Confidential & Proprietary

Company Confidential & Proprietary

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No. : 5  
Effective Date: 9/10/2012  
Page 143 of 165

- 2 - Method 8082, a five-point calibration is only analyzed for Aroclors 1016 and 1260.  
3 - This is a summary of the acceptance criteria, refer to the method SOP for specific or more information.  
4 - All abnormalities must be noted in an NCM.  
5 - Report all target compounds identified in the method blank above the MDL.  
6 - If unable to re-extract the samples because of insufficient sample volume or holding time has expired, then place a comment on the benchsheet and in LIMS.

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No. : 5  
Effective Date: 9/10/2012  
Page 144 of 165

Appendix 4: Summary of Calibration, QC and Corrective Action Procedures for GC Organics

Method	QC Check	Frequency	Acceptance Criteria <sup>3</sup>	Corrective Action <sup>4</sup>
EPA608	Minimum three-point (preferably five) initial calibration for all target analytes	Initial calibration prior to sample analysis. Perform instrument re-calibration once per year minimum.	RSD of CF $\leq 10\%$ Linear regression - correlation coefficient $r \geq 0.995$ , $r^2 \geq 0.990$ .	Correct problem then repeat initial calibration
	Initial calibration verification (ICV) must be from a 2 <sup>nd</sup> source	Immediately following initial calibration	All analytes within 15% of expected value	Correct problem then repeat initial calibration
	Continuing calibration verification (CCV)	Before sample analysis, and at the end of the analysis sequence	All analytes within <b>15% of expected value</b> and within the RTW. 608 (not 608.2*) must be within 15% of the true value.	Correct problem then repeat initial CCV (re-calibrate if necessary) and re-analyze all samples since last successful CCV.
	Breakdown check (Endrin and DDT)	Before sample analysis	Degradation $\leq 15\%$ for either Endrin or DDT.	Inlet/column maintenance; repeat breakdown check and re-analyze all samples since last successful breakdown check
	Method blank	One per analytical prep batch, <b>not to exceed 10 samples in a batch</b> .	No analytes detected $\geq \frac{1}{2}$ RL or MDL, whichever is greater <sup>5</sup>	Correct problem then re-prepare and analyze method blank and all samples processed with the contaminated blank
	LCS (QC check standard) must be from a 2 <sup>nd</sup> source.	One per prep batch, <b>not to exceed 10 samples in a batch</b> .	See Control Limits	Re-prepare and analyze the LCS and all samples in the affected analytical batch
	Surrogate(s)	Every sample, spiked sample, standard, and method blank	See Control Limits	Check system, re-inject, re-extract <sup>7</sup>
	MS must be from a 2 <sup>nd</sup> source.	One per batch per matrix, 10%, if insufficient sample for MS, then an additional LCS will be analyzed.	See Control Limits	All target compounds should be reported, and any compounds that are outside criteria must be within criteria in the LCS.
	Second-column confirmation	100% for all positive results	Same as for initial or primary column analysis	Same as for initial or primary column analysis. If the relative % difference of results between the 2 columns is greater than 40%, a comment should be placed in LIMS.
	Retention time window calculated for each analyte (see section 9 for how to calculate RTWs).	System set-up, with each new column or major instrument maintenance. Update the mid-RTW at the start of the run or as needed.	Each analyte of the LCS, MS/MSD and CCV must be within the calculated RTW.	Correct the problem and re-process or re-analyze samples. If questions, see the supervisor or Technical Manager.
	MDL verification	Minimum yearly	Detectable	Re-evaluate MDL standard used and MDL; see Technical Manager.

- 3 - This is a summary of the acceptance criteria, refer to the method SOP for specific or more information.  
4 - All abnormalities must be noted in an NCM.  
5 - 608.2 compounds are Methoxychlor and Permethrin.  
6 - Report all target compounds identified in the method blank above the MDL.  
7 - If unable to re-extract the samples because of insufficient sample volume or holding time has expired, then place a comment on the benchsheet and in LIMS.

Company Confidential & Proprietary

Company Confidential & Proprietary

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 5  
Effective Date: 9/10/2012  
Page 145 of 165

**Appendix 4: Summary of Calibration, QC and Corrective Action Procedures for GC/MS Organics**

Method	QC Check	Frequency	Acceptance Criteria <sup>a</sup>	Corrective Action <sup>a</sup>
SW8260B SW8270C SW8151mod	Check of mass spectral ion intensities <sup>b</sup> , i.e., Tune	Prior to initial calibration or Continuing calibration verification, every 12 hours	Refer to criteria listed in the method SOP for Tune criteria, including DDT, Benzidine and Pentachlorophenol requirements for 8270.	Retune the instrument and verify (instrument maintenance may be needed).
SW8260B SW8270C	Minimum five-point initial calibration for all target analytes	Initial calibration prior to sample analysis. Perform instrument re-calibration once per year minimum.	VOCs: SPCCs average RF $\geq 0.30$ or 0.1 depending on the compound and %RSD for RFs for CCCs $\leq 30\%$ and all other target analytes %RSD for RF $\leq 15\%$ .	Correct problem then repeat initial calibration
			SVOCs: SPCCs average RF $\geq 0.050$ and %RSD for RFs for CCCs $\leq 30\%$ and all other target analytes %RSD for RF $\leq 15\%$ , option (if %RSD is $> 15\%$ )-linear regression $r^2 \geq 0.99$ , $r \geq 0.995$ .	Correct problem then repeat initial calibration
			If the calibration is not considered linear by either %RSD or linear regression, then correct the problem and re-calibrate.	Correct problem then repeat initial calibration
	Initial calibration verification (ICV) must be from a 2 <sup>nd</sup> source.	Immediately following five-point initial calibration	All CCCs within 20% of expected value. <b>For DoD:</b> All analytes within 20% of expected value. See limits in SOP for poor-performing and other non-standard HSL compounds.	Correct problem then repeat initial calibration
	Relative Retention time window	Each sample	Relative retention time (RRT) of the analyte within 0.06 RRT units of the RRT of the internal standard	Correct problem then reprocess or re-analyze all samples analyzed since the last retention time check
	Continuing calibration verification (CCV)	Daily, before sample analysis and every 12 hours of analysis time	VOCs: SPCCs average RF $\geq 0.30$ or 0.1 depending on the compound;  <b>For DoD:</b> SPCCs average RF $\geq 0.050$ ; CCCs: $\leq 20\%$ difference (when using RFs) or drift (when using least squares regression).  See limits in SOP for poor-performing and other non-standard HSL compounds.	Correct problem then repeat initial calibration and re-analyze all samples since last successful CCV.
			SVOCs: SPCCs average RF $\geq 0.050$ ; and	

Company Confidential & Proprietary

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 5  
Effective Date: 9/10/2012  
Page 146 of 165

Method	QC Check	Frequency	Acceptance Criteria <sup>a</sup>	Corrective Action <sup>a</sup>
	Continuing calibration check		CCCs: $\leq 20\%$ difference (when using RFs) or drift (when using least squares regression) for standard list analytes, $\leq 35\%$ for non-standard analytes, and $\leq 55\%$ for poor performing analytes. <b>For DoD:</b> SPCCs average RF $\geq 0.050$ ; CCCs: $\leq 20\%$ difference (when using RFs) or drift (when using least squares regression). See limits in SOP for poor-performing and other non-standard HSL compounds.	
	Method blank	One per analytical prep batch	No analytes detected $\geq \frac{1}{2}$ RL or MDL, whichever is greater <sup>b</sup>	Correct problem then re-prep <sup>c</sup> and analyze method blank and all samples processed with the contaminated blank.
	Internal Standards	Every sample/standard and blank	Retention time $\sim 30$ seconds from retention time of the mid-point std. in the CCV/ICAL (sample/standard). EICP area within -50% to +100% of ICAL mid-point std for the CCV and -50% to +100% of the prior CCV for the samples. <b>For DoD:</b> EICP area within -50% to +100% of ICAL mid-point std for the CCV and samples.	Inspect mass spectrometer and GC for malfunctions; mandatory re-analysis of samples analyzed while system was malfunctioning (dilution of the sample may be required, see the supervisor or the Technical Manager for advice).
	LCS for all analytes must be from a 2 <sup>nd</sup> source.	One per prep batch, not to exceed the 20 samples in a batch.	See Control Limits Manual	Correct problem then re-prep <sup>c</sup> and analyze the LCS and all samples in the affected analytical batch.
	MS/MSD must be from a 2 <sup>nd</sup> source.	One per batch per matrix, if insufficient sample for MS/MSD, then a LCS/LCSD will be analyzed.	See Control Limits Manual	None (the LCS is used to evaluate to determine if the batch is acceptable).
	Surrogate(s)	Every sample, spike, standard, and blank	See Control Limits Manual	Check system, re-analyze, re-prep <sup>c</sup>
	MDL verification	Minimum yearly <b>For DOD:</b> Minimum quarterly	Detectable	Re-evaluate MDL standard used and MDL; see Technical Manager.
	LOQ verification	Minimum quarterly	Detectable	Re-evaluate LOQ standard used and LOQ; see Technical Manager.
SW8260B	pH check	All 8260 water samples.	pH $\leq 2$ .	If the pH is $> 2$ , then comment the data, in the PIPE database, and LIMS.
	Residual chlorine check (North Carolina samples only)	Each sample.	Residual chlorine should be negative.	If the residual chlorine is positive, then comment in the data and LIMS.

Company Confidential & Proprietary

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 5  
Effective Date: 9/10/2012  
Page 147 of 165

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 5  
Effective Date: 9/10/2012  
Page 148 of 165

- 1 - SW8260B requires BFB; SW8270C requires DFTPP
- 2 - This is a summary of the acceptance criteria, refer to the method SOP for specific or more information.
- 3 - All abnormalities must be noted in an NCM.
- 4 - Report all target compounds identified in the method blank above the MDL.
- 5 - If unable to re-prepare samples because of insufficient sample volume or the holding time has expired, then note in an NCM.

**Appendix 4: Summary of Calibration, QC and Corrective Action Procedures for GC/MS Organics**

Method	QC Check	Frequency	Acceptance Criteria <sup>2</sup>	Corrective Action <sup>3</sup>
EPA624 EPA625	Check of mass spectral ion intensities (i.e. Tune)	Prior to initial calibration or Continuing calibration verification every 12 hours.	Refer to criteria listed in the method SOP for Tune requirements including DDT, Benzidine and Pentachlorophenol criteria for 625.	Retune instrument and verify instrument maintenance may be needed.
	Minimum three-point initial calibration for all target analytes.	Initial calibration prior to sample analysis. Perform instrument re-calibration once per year minimum.	%RSD < 35%, if %RSD is > 35% then linear regression is used (for linear regression $r^2 \geq 0.99$ , $r \geq 0.995$ ).	If the calibration is not considered linear by either %RSD or linear regression, then correct problem then repeat initial calibration.
	Initial calibration verification (ICV), 20 ug/L, must be from a 2 <sup>nd</sup> source. May be the same as the LCS.	Immediately following initial calibration	See Control Limits Manual.	Correct problem then repeat initial calibration
	Relative Retention time window	Each sample	Retention time (RT) of the analyte within 30 seconds of the RT ( $\pm 0.25$ min. RTW is used) of the target.	Correct problem then reprocess or re-analyze all samples analyzed since the last retention time check.
EPA625	Continuing calibration verification (CCV)	Daily, before sample analysis and every 12 hours of analysis time.	All calibration analytes within 20% of expected value	Correct problem then repeat initial calibration and re-analyze all samples since last successful CCV.
EPA624 EPA625	Method blank	One per prep batch (not to exceed 20 samples per batch) or daily	See limits in SOP	Correct problem then re-prepare <sup>4</sup> and analyze method blank and all samples processed with the contaminated blank
	LCS for all analytes, 20 ug/L, must be from a 2 <sup>nd</sup> source. May be the same as the ICV.	One per prep batch (not to exceed 20 samples per batch) or daily	See Control Limits Manual	Correct problem then re-prepare <sup>4</sup> and analyze the LCS and all samples in the affected analytical batch.
	MS must be from a 2 <sup>nd</sup> source.	One per batch of 20 per matrix, if insufficient sample for MS, then a duplicate LCS will be analyzed.	See Control Limits Manual	All target compounds should be reported, and any compound that is outside criteria must be within criteria in the LCS.
	Surrogate(s)	Every sample, spiked sample, standard, and method blank	See Control Limits Manual	Correct problem then re-prepare <sup>4</sup> and analyze sample
EPA624 EPA625	Internal Standards	Every sample/standard	Retention time $\pm 30$ seconds from retention time of the mid-point std. in the CCV/ICAL (sample/standard). EICP area within -50% to +100% of ICAL mid-point std for the CCV and -50% to +100% of the prior CCV for the samples.	Inspect mass spectrometer and GC for malfunctions; mandatory re-analysis of samples analyzed while system was malfunctioning (dilution of the sample may be required, see the supervisor or the Technical Manager for advice).
EPA624	pH check	All 624 samples after analysis	pH should be $\leq 2$	If the pH is > 2, then comment the data, in the PIPE database, and LIMS.

Company Confidential & Proprietary

Company Confidential & Proprietary



Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 5  
Effective Date: 9/10/2012  
Page 149 of 165

Method	QC Check	Frequency	Acceptance Criteria <sup>1</sup>	Corrective Action <sup>2</sup>
EPA624	Residual chlorine check (North Carolina samples only)	All samples after analysis	Residual chlorine should be negative.	If the residual chlorine is positive, then comment the data and LIMS.
	MDL verification	Minimum yearly	Detectable	Re-evaluate MDL standard used and MDL; see Technical Manager.

- 1 - 624 requires BFB; 625 requires DFTPP  
2 - This is summary of the acceptance criteria, refer to the method SOP for specific or more information.  
3 - All abnormalities must be noted in an NCM.  
4 - Report all target compounds identified in the method blank above the MDL.  
5 - If unable to re-prepare samples because of insufficient sample volume or holding time has expired, then place a comment on the benchsheet and in LIMS.

Company Confidential & Proprietary

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 5  
Effective Date: 9/10/2012  
Page 150 of 165

Appendix 4: Summary of Calibration, QC and Corrective Action Procedures for GC Organics

Method	QC Check	Frequency	Acceptance Criteria <sup>1</sup>	Corrective Action <sup>2</sup>
SW8015 <sup>3</sup> SW8021	Five-point initial calibration for all target analytes	Initial calibration prior to sample analysis. Perform instrument re-calibration once per year minimum.	RSD of CF ≤ 20% Linear – least squares regression $r^2$ ≥ 0.99, $r$ ≥ 0.995	Correct problem then repeat initial calibration
	Initial calibration verification (ICV). must be from a 2 <sup>nd</sup> source.	Immediately following five-point initial calibration	All analytes within 15% of expected value	Correct problem then repeat initial calibration
	LCS for all analytes must be from a 2 <sup>nd</sup> source.	One per prep batch, not to exceed 20 samples in a batch.	See Control Limits	Re-prepare <sup>4</sup> and analyze the LCS and all samples in the affected analytical batch
	Continuing calibration verification (CCV)	Before sample analysis, after every 10 samples, and at the end of the analysis sequence	All analytes within 15% of expected value and within the RTW. <b>For DoD:</b> All analytes within 20% of expected value and within the RTW.	Correct problem then repeat initial CCV (re-calibrate if necessary) and re-analyze all samples since last successful CCV.
	Method blank	One per analytical prep batch, not to exceed 20 samples in a batch.	No analytes detected ≥ 1/2 RL or MDL, whichever is greater <sup>3</sup>	Correct problem then re-prepare <sup>4</sup> and analyze method blank and all samples processed with the contaminated blank
	Surrogate	Every sample, spiked sample, standard, and method blank	See Control Limits	Check system, re-analyze, re-prepare <sup>4</sup>
	MS/MSD must be from a 2 <sup>nd</sup> source.	One per batch per matrix, if insufficient sample for MS/MSD, then a LCS/LCSD will be analyzed.	See Control Limits	None (LCS is used to determine if data is acceptable).
	GC/MS confirmation	At the clients request or analyst judgment.		
	Retention time window calculated for each analyte (see section 9 for how to calculate RTWs).	System set-up, with each new column or major instrument maintenance. Update the mid-RTW as the start of the run or daily.	Each analyte of the LCS, MS/MSD and CCV must be within the calculated RTW.	Correct the problem and re-process or re-analyze samples. For questions, see the supervisor or Technical Manager.
	MDL verification	Minimum yearly <b>For DoD:</b> Minimum quarterly	Detectable	Re-evaluate MDL standard used and MDL; see Technical Manager.
	LOQ verification	<b>For DoD:</b> Minimum quarterly	Detectable	Re-evaluate LOQ standard used and LOQ; see Technical Manager.
SW8021	pH Check	All water samples after analysis.	pH should be less than 2.	If pH is > 2, then place a comment on the benchsheet and in LIMS.
MT-VPH MT-EPH WA-VPH WA-EPH NWTPH-EPH NWTPH-VPH	Five-point initial calibration for all target analytes	Initial calibration prior to sample analysis. Perform instrument re-calibration once per year minimum.	RSD of CF ≤ 20% Linear – least squares regression $r^2$ ≥ 0.99, $r$ ≥ 0.995	Correct problem then repeat initial calibration
	Initial calibration verification (ICV) or Petroleum Performance Check Standard	Immediately following five-point initial calibration	Must be within 70-130% recovery.	Correct problem then repeat initial calibration
	Retention time window calculated for each analyte (see section 9 for how to calculate RTWs).	System set-up, with each new column or major instrument maintenance. Update the mid-RTW as the start of the run or daily.	Each analyte of the LCS, MS/MSD and CCV must be within the calculated RTW.	Correct the problem and re-process or re-analyze samples. For questions, see the supervisor or Technical Manager.

Company Confidential & Proprietary

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 5  
Effective Date: 9/10/2012  
Page 151 of 165

Method	QC Check	Frequency	Acceptance Criteria <sup>1</sup>	Corrective Action <sup>2</sup>
	Continuing calibration verification (CCV)	Before sample analysis, after every 10 samples, and at the end of the analysis sequence	All analytes within 20% of expected value and within the RTW.	Correct problem then repeat initial CCV (re-calibrate if necessary) and re-analyze all samples since last successful CCV.
	Method blank	One per analytical prep batch, not to exceed 20 samples in a batch.	No analytes detected $\geq$ RL	Correct problem then re-prep <sup>3</sup> and analyze method blank and all samples processed with the contaminated blank
	LCS for all analytes must be from a 2 <sup>nd</sup> source	One per prep batch, not to exceed 20 samples in a batch.	Must be within 70-130% recovery.	Re-prep <sup>3</sup> and analyze the LCS and all samples in the affected analytical batch
	MS/MSD must be from a 2 <sup>nd</sup> source.	One per batch per matrix, if insufficient sample for MS/MSD, then a LCS/LCSD will be analyzed.	Must be within 50-150% recovery. For WA-VPH and WA-EPH: Must be within 70-130% recovery.	None (LCS is used to determine if data is acceptable).
	Duplicate	Each batch, less than 20	$\pm 25\%$	None
	Surrogate	Every sample, spiked sample, standard, and method blank.	Low end aliphatic compounds must be within 30-150% recovery and high end aliphatic compounds must be within 50-150% recovery. For WA-VPH and WA-EPH: Must be within 60-140% recovery.	Check system, re-analyze, re-prep <sup>3</sup>
	MDL and MDL verification	Minimum yearly	Detectable	Re-evaluate MDL standard used and MDL; see Technical Manager.
SW8011 EPA504.1	Five-point initial calibration for all target analytes (calibration standards should be prepped as the samples). Three point minimum for 504.1.	Initial calibration prior to sample analysis. Perform instrument re-calibration once per year minimum.	RSD of CF of each analyte $\leq 20\%$ <b>RSD of CF <math>&lt; 10\%</math> for Method 8011</b> Linear - $r^2 \geq 0.990$ , $r \geq 0.995$	Correct problem then repeat initial calibration
	Initial calibration verification (ICV) must be from a 2 <sup>nd</sup> source.	Immediately following initial calibration	All analytes within 20% of expected value	Correct problem then repeat initial calibration
	Continuing calibration verification (CCV)	Before sample analysis, after every 10 samples, and at the end of the analysis sequence	All analytes within 10% of expected value and within the RTW.	Repeat once. If subsequent CCV fails, correct problem, re-calibrate and re-analyze all samples since last successful CCV.
	Continuing calibration blank (CCB)	After each CCV	No analytes $\geq \frac{1}{2}$ RL	Correct problem and re-analyze all samples since last successful CCB, unless contamination is $< 10\times$ the conc present in the sample
	Method blank	One per analytical prep batch, not to exceed 10 samples in a batch.	No analytes $\geq \frac{1}{2}$ RL	Correct problem then re-prep <sup>3</sup> and analyze method blank and all samples processed with the contaminated blank
	LCS for all analytes must be from a 2 <sup>nd</sup> source.	One per prep batch, not to exceed 10 samples in a batch.	All analytes within 30% of expected value.	Re-prep <sup>3</sup> and analyze the LCS and all samples in the affected analytical batch

Company Confidential & Proprietary

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 5  
Effective Date: 9/10/2012  
Page 152 of 165

Method	QC Check	Frequency	Acceptance Criteria <sup>1</sup>	Corrective Action <sup>2</sup>
	Surrogate	Every sample, spike, standard, and method blank	See Control Limits	Check system, re-inject, re-extract <sup>4</sup>
	MS/MSD must be from a 2 <sup>nd</sup> source.	One per batch per matrix, if insufficient sample for MS/MSD, then a LCS/LCSD will be analyzed.	See Control Limits	(LCS is used to determine if data is acceptable). Add comments in LIMS.
	Second-column confirmation	100% for all positive results	Same as for initial or primary column analysis	Same as for initial or primary column analysis. If the relative % difference of results between the 2 columns is greater than 40%, a comment should be placed in LIMS.
	Retention time window calculated for each analyte (see section 20.1 for how to calculate RTW's).	System set-up, with each new column or major instrument maintenance. Update the mid-RTW as the start of the run or daily.	Each analyte of the LCS, MS/MSD and CCV must be within the calculated RTW.	Correct the problem and re-process or re-analyze samples. For questions, see the supervisor or Technical Manager.
	MDL / MDL verification	Minimum yearly For DOD: Minimum quarterly	MDL < Spike < 10X MDL Detected	Re-evaluate MDL standard used and MDL; see Technical Manager.
	LOQ verification	For DOD: Minimum quarterly	Detectable	Re-evaluate LOQ standard used and LOQ; see Technical Manager.

1 - This is a summary of the acceptance criteria, refer to the method SOP for specific or more information.

2 - All abnormalities must be noted on the data and in an NCM.

3 - Report all target compounds identified in the method blank above the MDL.

4 - If unable to re-prepare the samples because of insufficient sample volume or holding time has expired, then note in an NCM.

5 - For GRO and DRO, see state specific SOP/Method for acceptance criteria. If there is not a specific method for that state, then follow the acceptance criteria in this table.

Company Confidential & Proprietary

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 5  
Effective Date: 9/10/2012  
Page 153 of 165

**Appendix 4: Summary of Calibration, QC and Corrective Action Procedures for Method SW60108**

Method	QC Check	Frequency	Acceptance Criteria	Corrective Action
SW60108	Initial calibration (minimum 1 standard and a blank)	Daily initial calibration prior to sample analysis	N/A	N/A
	Second-source calibration verification (ICV)	Daily after initial calibration	All analytes within 10% of expected value	Correct problem then repeat initial calibration
	RL Standard (CRI)	Daily after initial calibration	All others: within 50% of expected value. For DoD: within 20% of expected value	Re-analyzed RL Standard, correct problem no samples may be analyzed without a valid RL Standard.
	Calibration blank (CB)	After every continuing calibration verification	Must be <RL. For DoD: Must be < LOD	Correct problem then analyze calibration blank and previous 10 samples
	Continuing calibration verification (CCV)	Before sample analysis, after every 10 samples, and at the end of the analysis sequence	All analytes within 10% of expected value and RSD of replicate integrations <5%	Repeat calibration and re-analyze all samples since last successful calibration
	Method blank	One per prep batch	No analytes detected $\geq \frac{1}{2}$ RL or MDL, whichever is greater <sup>1</sup> For DoD: Must be < $\frac{1}{2}$ RL	Correct problem then re-prepare and analyze method blank and all samples processed with the contaminated blank
	Interference check solution (ICS)	At the beginning of an analytical run	ICS-A: Absolute value of concentration for all non-spiked analytes < LOD (unless they are a verified trace impurity from one of the spike analytes) ICS-AB: Within 20% of expected value	Terminate analysis; correct problem; re-analyze ICS; re-analyze all affected samples
	LCS	One per prep batch	See Control Limits	Correct problem then re-prepare and analyze the LCS and all samples in the affected analytical batch
	MS/MSD	One per batch per matrix	See Control Limits	None
	Dilution test	Each new sample matrix	1:5 dilution must agree within 10% of the original determination	Perform post digestion spike addition
	Post digestion spike addition	When dilution test fails	Recovery within 25% of expected results	None
	MDL verification	Minimum yearly For DOD: Minimum quarterly	Detectable	Re-evaluate MDL standard used and MDL, see Technical Manager.
	LOQ verification	For DOD: Minimum quarterly	Detectable	Re-evaluate LOQ standard used and LOQ, see Technical Manager.

<sup>1</sup> - Report all targets identified in the method blank above the MDL.

Company Confidential & Proprietary

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 5  
Effective Date: 9/10/2012  
Page 154 of 165

**Appendix 4: Summary of Calibration, QC and Corrective Action Procedures for Method SW6020**

Method	QC Check	Frequency	Acceptance Criteria	Corrective Action
SW6020	Initial calibration (minimum 1 standard and a blank)	Daily initial calibration prior to sample analysis	If more than one calibration standard used, $r \geq 0.995$	Correct problem then repeat initial calibration
	Second-source calibration verification (ICV)	Daily after initial calibration	All analytes within 10% of expected value	Correct problem then repeat initial calibration
	RL Standard (CRI)	Daily after initial calibration	All others: within 50% of expected value. For DoD: within 20% of expected value	Re-analyzed RL Standard, correct problem no samples may be analyzed without a valid RL Standard.
	Calibration blank (CB)	After every continuing calibration verification	Must be <RL. For DoD: Must be < LOD	Correct problem then analyze calibration blank and previous 10 samples
	Continuing calibration verification (CCV)	Before sample analysis, after every 10 samples, and at the end of the analysis sequence	All analytes within 10% of expected value and RSD of replicate integrations <5%	Repeat calibration and re-analyze all samples since last successful calibration
	Linear Dynamic Range	Every 6 months	Within 10% of expected value	N/A
	Method blank	One per prep batch	No analytes detected $\geq$ RL or MDL, whichever is greater <sup>1</sup> For DoD: Must be < $\frac{1}{2}$ RL	Correct problem then re-prepare and analyze method blank and all samples processed with the contaminated blank
	Interference check solution (ICS)	At the beginning of an analytical run	ICS-A: Absolute value of concentration for all non-spiked analytes < LOD (unless they are a verified trace impurity from one of the spike analytes) ICS-AB: Within 20% of expected value	Terminate analysis; correct problem; re-analyze ICS; re-analyze all affected samples
	Internal standards (IS)	Every sample	IS intensity within 30-120% of intensity of the IS in the initial calibration	Perform corrective action as described in 6020.
	RL standard	Daily, after one point initial calibration	Within 20% of expected value	Correct problem then re-analyze samples
	LCS	One per prep batch	See Control Limits	Correct problem then re-prepare and analyze the LCS and all samples in the affected analytical batch
	MS/MSD	One per batch per matrix	See Control Limits	None
	Dilution test	Each new sample matrix	1:5 dilution must agree within 10% of the original determination	Perform post digestion spike addition

Post digestion spike addition	When dilution test fails	Recovery within 25% of expected results	None.
MDL verification	Minimum yearly For DOD: Minimum quarterly	Detectable	Re-evaluate MDL standard used and MDL, see Technical Manager.
LOQ verification	For DOD: Minimum quarterly	Detectable	Re-evaluate LOQ standard used and LOQ, see Technical Manager.

Company Confidential & Proprietary

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 5  
Effective Date: 9/10/2012  
Page 155 of 165

1 - Report all targets identified in the method blank above the MDL.

**Appendix 4: Summary of Calibration, QC and Corrective Action Procedures for Method SW7196**

Method	QC Check	Frequency	Acceptance Criteria	Corrective Action
SW7196A	Initial calibration (minimum three standards and a blank)	Initial calibration prior to sample analysis	$r^2 \geq 0.99$ , $r > 0.995$ for linear regression	Correct problem then repeat initial calibration
	Second-source calibration verification (ICV)	Immediately following initial calibration	All analytes within 10% of expected value	Correct problem then repeat initial calibration
	Continuing calibration verification (CCV)	Beginning and after every 10 samples and at the end of the analysis sequence	All analytes within 20% of expected value	Correct problem then repeat initial calibration and re-analyze all samples since last successful calibration
	Verification check to ensure lack of reducing condition and/or interference	Once for every sample matrix analyzed	Spike recovery between 85-115%	If check indicates interference, dilute and re-analyze sample persistent interference indicates the need to use and alternate method
	Method blank	One per prep batch	No analytes detected $\geq \frac{1}{2}$ RL or MDL, whichever is greater <sup>1</sup>	Correct problem then re-prepare and analyze method blank and all samples processed with the contaminated blank
	MS/MSD	One per 20 samples per matrix	See Control Limits Manual	none
	LCS	One per batch	See Control Limits Manual	Re-prepare, re-analyze all affected samples
	MDL verification	Minimum yearly	Detectable	Re-evaluate MDL standard used and MDL; see Technical Manager

1 - Report all targets identified in the method blank above the MDL.

Company Confidential & Proprietary

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 5  
Effective Date: 9/10/2012  
Page 156 of 165

**Appendix 4: Summary of Calibration, QC and Corrective Action Procedures for Method SW7470/SW7471**

Method	QC Check	Frequency	Acceptance Criteria	Corrective Action
SW7470A SW7471A	Initial calibration (minimum 5 standards and a blank)	Daily initial calibration prior to sample analysis. Perform instrument re-calibration once per year minimum	$r^2 \geq 0.99$ , $r \geq 0.995$ for linear regression	Correct problem then repeat initial calibration. If calibration fails again, re-digest the entire digestion batch.
	Second-source calibration verification (ICV)	Immediately following initial daily calibration	Analytes within 10% of expected value	Correct problem then repeat initial calibration. If calibration fails again, re-digest the entire digestion batch.
	Calibration blank	Once per initial daily calibration	No analytes detected > MDL	Correct problem then re-digest and re-analyze calibration and entire digestion batch
	Continuing calibration verification (CCV)	Before sample analysis, after every 10 samples, and at the end of the analysis sequence	Analytes within 20% of expected value	Correct problem then repeat all QC and samples since last successful calibration. If the CCV fails again upon reanalysis, reprep the entire digestion batch.
	Method blank	One per prep batch	No analytes detected > RL <b>For DoD:</b> Must be < $\frac{1}{2}$ RL	Correct problem then re-prepare and analyze method blank, all samples, and QC processed with the contaminated blank
	LCS	One per prep batch	See Control Limits Manual	Correct problem then re-prepare and analyze the LCS, all samples, and QC in the affected analytical batch
	MS/MSD	One per batch per matrix	See Control Limits Manual	None
	MDL verification	Minimum yearly <b>For DOD:</b> Minimum quarterly	Detectable	Re-evaluate MDL standard used and MDL; see Technical Manager
	LOQ verification	<b>For DOD:</b> Minimum quarterly	Detectable	Re-evaluate LOQ standard used and LOQ; see Technical Manager

1 - Report all targets identified in the method blank above the MDL.

Company Confidential & Proprietary

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 5  
Effective Date: 9/10/2012  
Page 157 of 165

**Appendix 4: Summary of Calibration, QC and Corrective Action Procedures for Method SW9010/SW9012/SW9014**

Method	QC Check	Frequency	Acceptance Criteria	Corrective Action
SW9010 SW9012A SW9014	Initial calibration (six standards and a calibration blank)	Initial daily calibration prior to sample analysis. Perform instrument re-calibration once per year minimum.	$r^2 \geq 0.99$ , $r \geq 0.995$ for linear regression	Correct problem then repeat initial calibration
	Distilled standards (one high and one low)	Once per calibration	Analyte within 10% of true value. <b>For DoD:</b> Analyte within 15% of true value.	Correct problem then repeat distilled standards
	Second-source calibration verification (ICV)	Immediately following initial daily calibration	Analyte within 15% of expected value.	Correct problem then repeat initial calibration
	Continuing calibration verification (CCV)	Beginning and after every 10 samples and at the end of the analysis sequence	Analyte within 15% of expected value	Correct problem then repeat initial Continuing calibration verification and re-analyze all samples since last successful Continuing calibration verification.
	Calibration blank	Once per initial daily calibration	No analytes detected $\geq$ RL <b>For DoD:</b> Must be $< \frac{1}{2}$ RL	Correct problem then re-analyze calibration blank and all samples associated with blank
	Method blank	One per prep batch	No analyte detected $\geq$ RL or MDL, whichever is greater <b>For DoD:</b> Must be $< \frac{1}{2}$ RL	Correct problem then re-prepare and analyze method blank and all samples processed with the contaminated blank
	LCS	One per batch per matrix	See Control Limits in LIMS	Re-prepare, re-run affected samples
	MS/MSD	One per batch per matrix	See Control Limits in LIMS	None
	MDL verification	Minimum quarterly <b>For DOD:</b> Minimum quarterly	Detectable	Re-evaluate MDL standard used and MDL; see Technical Manager.
	LOQ verification	<b>For DOD:</b> Minimum quarterly	Detectable	Re-evaluate LOQ standard used and LOQ; see Technical Manager.

Company Confidential & Proprietary

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 5  
Effective Date: 9/10/2012  
Page 158 of 165

**Appendix 4: Summary of Calibration, QC and Corrective Action Procedures for Mercury**

Method	QC Check	Frequency	Acceptance Criteria	Corrective Action
EPA245.1	Initial calibration (minimum 5 standards and a blank)	Daily initial calibration prior to sample analysis. Perform instrument re-calibration once per year minimum.	$r^2 \geq 0.99$ , $r \geq 0.995$ for linear regression	Correct problem then repeat initial calibration
	Second-source calibration verification (ICV)	Immediately following five-point initial calibration	Analyte within 5% of expected value	Correct problem then repeat initial calibration
	Calibration blank	Once per initial daily calibration	No analytes detected $\geq$ MDL	Correct problem then re-analyze calibration blank and all samples associated with blank
	Continuing calibration verification (CCV)	Before sample analysis, after every 10 samples, and at the end of the analysis sequence	Analyte within 10% of true value	Correct problem then repeat calibration and re-analyze all samples and QC since last successful calibration
	LCS	One per prep batch	All analytes within 15% of expected value	Correct problem then re-prepare and analyze the LCS, all samples, and QC in the affected analytical batch
	Matrix Spike/Matrix Spike Duplicate	One per batch or 10 samples	All analytes within 25% of expected value	None
	Method Blank	One per batch	No analytes $> \frac{1}{2}$ RL	Reprep
	MDL verification	Minimum yearly	Detectable	Re-evaluate MDL standard used and MDL; see Technical Manager.

1 - Report all targets identified in the method blank above the MDL.

Company Confidential & Proprietary

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 5  
Effective Date: 9/10/2012  
Page 159 of 165

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 5  
Effective Date: 9/10/2012  
Page 160 of 165

Appendix 4: Summary of Calibration, QC and Corrective Action Procedures for ICP Metals

Method	QC Check	Frequency	Acceptance Criteria	Corrective Action
EPA200.7	Initial calibration (minimum 1 standard and a blank)	Daily initial calibration prior to sample analysis.	N/A	N/A
	Second-source calibration verification (ICV/QCS)	Each calibration	Value of all analytes within 5% of expected value	Correct problem then repeat initial calibration
	Linear Dynamic Range	Once annually	All analytes within 10% of expected value	Calibration range lowered to meet LDR results
	Calibration blank	After every Continuing calibration verification	No analytes detected (< IDL or MDL, whichever is greater)	Correct problem then analyze calibration blank and previous 10 samples
	Continuing calibration verification (CCV/IPC)	Before sample analysis, after every 10 samples, and at the end of the analysis sequence	All analytes within 10%	Repeat calibration and re-analyze all samples since last successful calibration
	Method blank	One per prep batch	No analytes detected, < 10% of sample conc OR < 2.2 x MDL <sup>1</sup>	Correct problem then re-prepare and analyze method blank and all samples processed with the contaminated blank
	Interference check solution (ICSA/ACSAB)	At the beginning of an analytical run, daily	Monitor for interferences.	Terminate analysis; correct problem; re-analyze ICS; re-analyze all affected samples
	LCS/LFB	One per prep batch	All analytes within 15% of expected value	Correct problem then re-prepare and analyze the LCS and all samples in the affected analytical batch
	Dilution test	Each new sample matrix	1:5 dilution must agree within 10% of the original determination for concentrations > 50 X IDL	Perform post digestion spike addition
	Post digestion spike addition	When dilution test fails	Recovery within 15% of expected results	Correct problem then re-analyze post digestion spike addition
	Matrix Spike/Matrix Spike Duplicate	One per batch of 10 samples	All analytes within 30% of expected value	None
	MDL verification	Minimum yearly	Detectable	Re-evaluate MDL standard used and MDL; see Technical Manager.

1 - Report all targets identified in the method blank above the MDL.

Company Confidential & Proprietary

Appendix 4: Summary of Calibration, QC and Corrective Action Procedures for ICP-MS Metals

Method	QC Check	Frequency	Acceptance Criteria	Corrective Action
EPA200.8	Initial calibration (minimum 1 standard and a blank)	Daily initial calibration prior to sample analysis.	N/A	N/A
	Second-source calibration verification (ICV/QCS)	Each calibration	Value of all analytes within 10% of expected value	Correct problem then repeat initial calibration
	Calibration blank (ICB)	After every Continuing calibration verification	No analytes detected (< RL)	Correct problem then analyze calibration blank and previous 10 samples
	Continuing calibration verification (CCV)	Before sample analysis, after every 10 samples, and at the end of the analysis sequence	All analytes within 10%	Repeat calibration and re-analyze all samples since last successful calibration
	Method Blank/Laboratory Reagent Blank	One per prep batch	No analytes detected, < 10% of sample conc. OR < 2.2xMDL <sup>1</sup>	Correct problem then re-prepare and analyze method blank and all samples processed with the contaminated blank
	Interference check solution (ICSA/ICSAB)	At the beginning of an analytical run, daily	Monitor for possible interferences.	Terminate analysis; correct problem; re-analyze ICS; re-analyze all affected samples
	LCS/LFB	One per prep batch	All analytes within 15% of expected value	Correct problem then re-prepare and analyze the LCS and all samples in the affected analytical batch
	Internal Standard	Every sample	IS intensity within 60-125% of intensity of the IS in the initial calibration	Perform corrective action as described in 200.8.
	Matrix Spike/Matrix Spike Duplicate (MS/MSD)	One per batch of 10 samples	All analytes within 30% of expected value	None
	MDL verification	Minimum yearly	Detectable	Re-evaluate MDL standard used and MDL; see Technical Manager.

1 - Report all targets identified in the method blank above the MDL.

Company Confidential & Proprietary

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 5  
Effective Date: 9/10/2012  
Page 161 of 165

**Appendix 4: Summary of Calibration, QC and Corrective Action Procedures for Gravimetric Analyses**

Method	QC Check	Frequency	Acceptance Criteria	Corrective Action
EPA160.1 (TDS)	Verification standard– single standard (if available)	Each batch	±10%	Repeat
SM2540 C (TDS)	Method blank	Each batch	No analytes detected ≥ RL	Repeat, rerun
EPA160.2 (TSS)				
SM2540D (TSS)				
EPA160.3 (TS)				
SM2540B (TS)				
EPA160.4 (TVS)				
	Duplicate	Each batch, less than 20	±20%	None

Company Confidential & Proprietary

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 5  
Effective Date: 9/10/2012  
Page 162 of 165

**Appendix 4: Summary of Calibration, QC and Corrective Action Procedures for Titrimetric Analyses**

Method	QC Check	Frequency	Acceptance Criteria	Corrective Action
EPA310.1: Alkalinity, SM2320: HCO <sub>3</sub> <sup>-</sup> , CO <sub>3</sub> <sup>2-</sup>	Verification standard– single standard (if available)	Each batch	±10%	Repeat, check
	Method blank	Each batch	No analyte detected ≥ RL	Repeat batch
	Duplicate	Each batch	±20%	None

**Appendix 4: Summary of Calibration, QC and Corrective Action Procedures for Electrometric Analyses**

Method	QC Check	Frequency	Acceptance Criteria	Corrective Action
EPA405.1: BOD <sup>1</sup> , CBOD <sup>1</sup>	Calibration Curve – minimum of 5 standards	Initial Calibration. Perform re-calibration once per year minimum	±10%, r <sup>2</sup> ≥ 0.99, r ≥ 0.995	Recalibrate
EPA120.1: Cond <sup>1</sup> , SM2510B: Cond <sup>1</sup>	Independent calibration verification (second source) (ICV)	Immediately after initial calibration	±10%	Recalibrate
EPA150.1: pH, SW9040B, 9045C.pH, EPA180.1: Turbidity	Continuing calibration verification (CCV)	Beginning, every 10 samples, and end of batch	±10%	Rerun
	Method blank	Each batch	No analyte detected ≥ ½ report limit or MDL, whichever is greater <sup>2</sup>	Reprep
	LCS	Each batch	±10%	Rerun batch
	MS/MSD	Each batch	± 20%	None
	Duplicate	When spike not available	±20%	None
	MDL verification	Minimum yearly	Detectable	Re-evaluate MDL standard used and MDL; see Technical Manager.

<sup>1</sup> - Calibration curve does not apply.

<sup>2</sup> - Report all targets identified in the method blank above the MDL.

Company Confidential & Proprietary



Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 5  
Effective Date: 9/10/2012  
Page 163 of 165

**Appendix 4: Summary of Calibration, QC and Corrective Action Procedures for Spectrophotometric Analyses**

Method	QC Check	Frequency	Acceptance Criteria	Corrective Action
EPA350.1: NH <sub>3</sub> EPA410.4: COD	Calibration curve – minimum 5 point	Initial. Perform re-calibration once per year minimum	RSD <10%, $r^2$ > 0.99, $r$ ≥ 0.995	Recalibrate
EPA335.4: CN EPA353.2 NO <sub>2</sub> /NO <sub>3</sub> EPA365.1 SM4500-P E T, Phos EPA415.1: TOC SM4500-CN1: WAD CH EPA9060: TOC SM4500-Cr D Hexchrome SM4500-N <sub>org</sub> C	Independent calibration verification – mid-level, second-source required (ICV)	Immediately following initial calibration.	See Control Limits in LIMS	Recalibrate
	Continuing calibration verification (CCV)	Beginning, every 10 samples, and at end of sequence	See Control Limits in LIMS	Correct, recalibrate
	Method blank	Each use	No analyte detected ≥ ½ report limit or MDL, whichever is greater <sup>1</sup>	Reprep, rerun
	MS/MSD	Each batch, less than 20	See Control Limits in LIMS	None
	LCS	Each batch	See Control Limits in LIMS	Rerun
	MDL verification	Minimum yearly	Detectable	Re-evaluate MDL standard used and MDL; see Technical Manager.

Company Confidential & Proprietary

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 5  
Effective Date: 9/10/2012  
Page 164 of 165

**Appendix 4: Summary of Calibration, QC and Corrective Action Procedures for Ion Chromatographic Analyses**

Method	QC Check	Frequency	Acceptance Criteria	Corrective Action
EPA300: SW9056A: Bromide Chloride Fluoride Nitrate Nitrite Sulfate.	Calibration Curve – Minimum 5-point calibration	Initial calibration. Perform instrument re-calibration once per year minimum.	RSD ± 10%, $r^2$ > 0.99, $r$ ≥ 0.995.	Recalibrate
	Calibration verification (ICV), second source	Immediately following initial calibration	±10%	Recalibrate
	Continuing calibration verification (CCV)	Each use, beginning, every 10 samples, end of batch	± 10%	Rerun affected samples
	Method blank	Each batch	No analyte detected ≥ ½ report limit or MDL, whichever is greater <sup>1</sup>	Rerun batch
	LCS	Each batch	±10%	Rerun batch
	MS/MSD <sup>1</sup>	Each batch	±20%	None, use LCS
	Duplicate <sup>2</sup>	Each batch	±30%	None
	MDL verification	Minimum yearly For DOD: Minimum quarterly	Detectable	Re-evaluate MDL standard used and MDL; see Technical Manager.
	LOQ verification	For DOD: Minimum quarterly	Detectable	Re-evaluate LOQ standard used and LOQ; see Technical Manager.

1 - Only applies to EPA300, SW9056.

2 - Oil only.

3 - Report all targets identified in the method blank above the MDL.

**Appendix 4: Summary of Calibration, QC and Corrective Action Procedures for Physical Analyses**

Method	QC Check	Frequency	Acceptance Criteria	Corrective Action
SW1020A: Flash Point EPA160.5: Settleable Solids	Method blank	Each batch (if applicable to the method)	No analyte detected ≥ ½ report limit or MDL, whichever is greater <sup>1</sup>	Repeat, rerun
	Two standards for Flash Point One standard for BTU 1 Known for Settleable Solids	Each batch	See Control Limits in LIMS	Rerun batch
	Duplicate	Each batch	See Control Limits in LIMS	None

1 - Report all targets identified in the method blank above the MDL.

Company Confidential & Proprietary

Document Uncontrolled When Printed

Document No. TA-QAM  
Revision No.: 5  
Effective Date: 9/10/2012  
Page 165 of 165

**Appendix 4: Summary of Calibration, QC and Corrective Action Procedures for Oil & Grease Analyses**

Method	QC Check	Frequency	Acceptance Criteria	Corrective Action
EPA1664A	Verification standard	Single standard	±10% PAR standard	Rerun
	Method blank	Each batch	No analyte detected ≥ report limit	Repeat batch
	LCS	Each batch	See Control Limits	Repeat batch
	Duplicate	Each batch	See Control Limits	
	MS/MSD	Each batch	See Control Limits	None, use LCS
	LOQ verification	For DOD: Minimum quarterly	Detectable	Re-evaluate LOQ standard used and LOQ; see Technical Manager.

Company Confidential & Proprietary

Analysis Group Description	Method Description	Method Code
Soil - Dalapon	Herbicides (GC/MS)	8151A_MS

Analyte Description	CAS Number	Reference RL - Limit
2,4-Dichlorophenylacetic acid	19719-28-9	
Dalapon	75-99-0	5.00

Soil - Dalapon	Extraction (Herbicides)	8151A_SP
----------------	-------------------------	----------

Reference RL - Units	Reference MDL - Limit	Reference MDL - Units	Reference LCSREC - Recovery Low	Reference LCSREC - Recovery High	Reference LCSREC - Units
ug/Kg		ug/Kg			%
ug/Kg	1.50	ug/Kg	23	117	%

Reference LCSRPD - Precision	Reference LCSRPD - Units	Reference MSREC - Recovery Low	Reference MSREC - Recovery High	Reference MSREC - Units
	%			%
30	%	23	117	%

Reference MSRPD - Precision	Reference MSRPD - Units	Reference SUREC - Recovery Low	Reference SUREC - Recovery High	Reference SUREC - Units
	%	51	129	%
30	%			%