
Microfiche No.					
	OTS0534	398 .			
New Doc I.D.		<u> </u>	Old Doc I.D.		
	86-920000193				
Date Produced		Date Recieve		TSCA section	
	11/01/80	Date Nacion	11/27/91	130A Section	BD .
Submitting Organiz	- HIGH				
	MOE	AY CORP	**		
Contractor					
	INVERESK RES	INTL LTD			
Document Title					
TESTING O	F THE CELL TRA	NSFORMAT	ION ACTIVITY OF	HE 1003	
	PORT) WITH COV				
	*				
Chemical Ca		d ENE TOO	DVANATE (DOLG 07	-0\	
POLYMETHY	LENE POLYPHENY	LENE 180	CYANATE (9016-87	-3)	

# CONTAINS NO CE

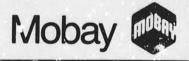

91 NOV 27 AM 9: 38

A Bayer USA INC COMPANY

Mobay Corporation Health, Environment, Safety & Plant Management

Mobay Road Pittsburgh, PA 15205-9741 Phone: 412 777-2000

November 22, 1931

Document Processing Center TS-790 Office of Toxic Substances Room L-100 Environmental Protection Agency 401 M Street SW Washington, DC 20460

Attention: 8(d) Health and Safety Reporting Rule

(Notification/Reporting)

#### Gentlemen:

Enclosed are copies of Health and Safety Studies submitted on behalf of Mobay Corporation, Mobay Road, Pittsburgh, Pennsylvania 15205. We are filing these Health and Safety Studies to comply with the regulations codified at 40 CFR, Part 716. This submission contains no Confidential Business Information (CBI).

The information required at 40 CFR 716.30 is given below.

Polymethylene polyphenylene isocyanate Chemical Name:

> CAS No: 9016-87-9

Testing of the Cell Transformation Activity of HE 1003 Name of Study:

Francis J. Rattay Submitting Official:

> Manager, Regulatory Compliance Title:

Mobay Road Address:

Pittsburgh, Pa 15205

(412) 777-7471 (412) 777-7484 Telephone No.:

FAX No.:

If you have any questions, please contact me.

86920000 193

Sincerely,

Francis J. Rettav

Manager, Regulatory Compliance

(412) 777-7471

Attachments

Certified Mail No.: P 276 377 332

# CONTAINS NO CBI

TESTING OF THE CELL TRANSFORMATION ACTIVITY OF HE 1003

IRI Project No. 703863

Bayer Study No. HE 1003/013



)

IRI Inveresk Research International

INVERESK RESEARCH INTERNATIONAL Report No. 1863

# CONFIDENTIAL

# TESTING OF THE CELL TRANSFORMATION ACTIVITY OF HE 1003

IRI Project No. 703863

Bayer Study No. HE 1003/013

Authors:

A. Poole W.J. Harris

To:

Issued by:

Bayer AG Inveresk
Institut fur Toxikologie Musselbur
Friedrich-Ebart-Strasse 217-319 EH21 7UB
Wuppertal 1 Scotland
West Germany

Inveresk Research International Musselburgh EH21 7UB Scotland

November 1980

Testing of the cell transformation activity of HE 1003

IRI Project No. 70 38 63

Bayer Study No. HE 1003/013

Bei HE 1003 handelt es sich um Desmodur 44 V 20. Desmodur 44 V 20 (MDI) ist ein technisches Gemisch mit der Summenformel  $^{\rm C}_{15}{}^{\rm H}_{10}{}^{\rm O}_{2}{}^{\rm N}_{2}$  und der Struktur

$$NCO - \langle \overline{O} \rangle - CH_2 - \langle \overline{O} \rangle - OCN$$

Die Untersuchung trug die BAYER Studie Nr.:

DESMODUR 44 V 20/013

Wuppertal, den 22. September 1981

Dr. B. Herbold

Monitor

#### AUTHENTICATION

"I, the undersigned, hereby declare that this work was performed under my supervision, according to the procedures herein described and that this report represents a true and accurate record of the results obtained."

W.J. Harris, B.Sc., Ph.D.
Principal Investigator

Com Expor

Project No. 703863

Report No. 1863

# QUALITY ASSURANCE AUTHENTICATION

The execution of this type of short-term study is not individually inspected. The processes involved are inspected at intervals according to a pre-determined schedule.

The report has been audited by IRI Quality Assurance personnel according to the appropriate Standard Operating Procedure and is considered to describe the methods and procedures used in the study. The reported results accurately reflect the original data of the study.

IRI PROJECT NO. 703863

Report No. 1863

Quality Assurance Manager

Date 3/4 august 1981.

# CONTENTS

	Page
STAFF INVOLVEMENT	1
SUMMARY	2
INTRODUCTION	3
MATERIALS AND METHODS	4
RESULTS AND DISCUSSION	11
CONCLUSION	12
REFERENCES	13
TABLES 1-8	14
FIGURES 1-7	22
DINAL DACE OF DEDODE	28

# STAFF INVOLVEMENT

Principal Investigator: W.J. Harris, B.Sc., Ph.D.

Project Leader: A. Poole, B.Sc., Ph.D.

Experimental Assistance: N. Hunter, H.N.D.

C. Green, B.Sc.

Quality Assurance: A. Waddell, B.Sc., Ph.D.

#### SUMMARY

The compound HE 1003 was tested for potential carcinogenicity in a cell culture transformation assay.

Cells treated with HE 1003, in the presence of S-9, caused at the  $LC_{50}$ , 47.2, 13.2, 1.8 and 10.2 fold increases in transformation frequency as compared to the negative controls. In the absence of S-9 3.1, 5.6 and 13.5 fold increases were observed.

The significant increases in the transformation frequencies at the LC<sub>50</sub> were also accompanied by increases in absolute numbers of transformed colonies in cells treated with HE 1003.

By the criterion used in this laboratory HE 1003 is considered to be a potential cell transforming agent.

#### INTRODUCTION

HE 1003 was tested for potential carcinogenicity by a cell culture transformation assay.

The assay is based on the observations that, following exposure to a carcinogen, malignantly transformed cells unlike their non-malignant counterparts, will undergo sufficient divisions in a soft agar medium to produce macroscopically visible colonies (Shin, et al 1975 and di Mayorca et al 1973).

As normal laboratory stocks of BHK 21 C13 cells contain spontaneously transformed variants which frequently cause a high background of macroscopic colonies when plated in soft agar, i.e.  $50/10^6$  plated cells, a subline selected from normal stocks able to express the transformed phenotype in response to treatment with chemical carcinogens, but with a low spontaneous transformation rate, was used in these experiments. The procedure used otherwise is a modification of that described by Styles (1977).

This work was carried out at the Inveresk Gate Laboratories of Inveresk Research International, Musselburgh, EH21 7UB between August and October 1980.

#### MATERIALS AND METHODS

Sterile procedures were used throughout preparation of materials and experimental methods.

#### Chemicals

The brown liquid labelled HE 1003 was received from Bayer AG, Wuppertal, West Germany on the 15 August 1980. The sample was stored in the dark at room temperature.

The positive control substances N-methyl-N'-nitro-N-nitrosoguanidine (MNNG) and 2-acetamidofluorene (2-AAF) were obtained from Sigma London Chemical Company Limited, U.K. and Noch-Light Laboratories, Colnbrook, Bucks., respectively.

The polychlorinated biphenyl mixture, Arcclor 1254, was received from Analabs Incorporated, Newhaven, Connecticut, U.S.A.

#### Cells

These were a subline from a sample of BHK 21 C13 of Syrian hamster kidney cells provided by the Imperial Cancer Research Fund.

#### Culture Media

For cell cultivation, the medium was Dulbecco's Modification of Eagle's Minimum Essential Medium (DMEM), with a concentration of sodium bicarbonate of 2 g/l to permit equilibration with a gas phase of 5%  $\rm CO_2$  in air, it was supplemented with newborn calf serum (10%, v/v), selected by pre-screening a number of serum samples for cloning efficiency, and contained gentamicin (50  $\mu$ g/ml).

For incubation with the compound under test, the cells, in order to prevent clumping during exposure, were suspended in Eagle's Minimum Essential Medium modified for suspension cultures (MEMS), buffered with 20 mM HEPES, the S-9 mixture was added to ascertain if the carcinogenic potential of the compound was affected by the metabolic activity of liver enzymes

Culture media and sera were obtained from Flow Laboratories, Irvine, Scotland and other chemicals from Sigma London Chemical Company Limited.

#### Maintenance of Stock Cultures

These were grown as monolayers in Nunc flasks and, in order to minimise selection for spontaneously transformed variants which accumulate when stocks reach confluence (Kao et al 1975), the cells were never grown to confluency.

For subcultivation, the medium was removed and the cells treated for one minute in a solution of 0.25% trypsin in phosphate buffered balanced salt solution containing EDTA (0.002% w/v). After removal of the enzyme solution, the flasks were left to incubate at 37°C until the cells began to detach from the plastic. 5 ml of fresh culture medium was then added and the cells brought into suspension by repeated aspiration through a sterile 10 ml pipette. Aliquots of the cell suspension were then added to medium in fresh culture flasks, the usual ratio for division of monolayers being 1:30.

#### Preparation of S-9 Mix

#### Animal Treatment

Male rats of the Bantin and Kingman Fischer strain weighing between 200 and 300 g were injected (i.p.) with Aroclor 1254

diluted in corn oil to a concentration of 200 mg/ml at a dosage of 500 mg/kg body weight to induce microsomal enzyme activity.

The animals were killed by cervical dislocation 5 days after treatment following a 16 h fasting period.

# Preparation of the 9,000 g Supernatant Fluid from Livers

Under aseptic conditions, livers from the freshly killed animals were carefully removed and weighed in sterile beakers containing ice-cold 0.15 M-KCl. Further ice-cold salt solution was added to the beakers to give a final volume equivalent to 3 times the weight of the livers which were subsequently finely chopped with sterile long-handled scissors before being transferred to a Potter homogeniser.

The chopped livers were homogenised in a sterile glass assel by 8 strokes of a Teflon pestle rotating at about 1,200 r.p.m. The homogenate was then transferred to sterile polypropylene tubes and centrifuged at 9,000 g for 10 min at 0°C.

The supernatant fraction was decanted into sterile containers and stored in liquid nitrogen until required.

#### 5-9

Under aseptic conditions, the S-9 mix was prepared as follows:

Time-cold 0.05 M-phosphate buffer, pH 7.4, was added to preweighed co-factors to give a final concentration in the S-9 mix of:

NADP-di-Na-salt 4 mM (= 3.366 mg/ml) Glucose-6-phosphate-di-Na-salt 5 mM (= 1.521 mg/ml)  $MgCl_2.6H_2O$  8 mM (= 1.626 mg/ml) KCl 33 mM (= 2.460 mg/ml)

The solution was sterilised by passage through a 0.22  $\mu m$  filter and mixed with the liver 9,000 g supernatant fluid in the following proportions:

co-factor 1 part liver preparation 1 part

# Preliminary Toxicity Test

This was done to establish the range of concentrations to be used in the cell transformation assay.

The cells were harvested and suspended in growth medium as for subculture, sedimented by centrifugation at 200 g for 5 min and resuspended in MEMS at a density of 10<sup>6</sup> cells/ml. 1 ml samples were then pipetted into plastic universal bottles (Sterilin Limited).

The test material was dissolved in dimethylsulphoxide at concentrations of 100, 10, 1.0, 0.1 and 0.01 mg/ml and 10  $\mu$ l samples were added to duplicate cell suspensions to give final concentrations of 1,000, 100, 10, 1.0 and 0.1  $\mu$ g/ml.

After incubation for 4 h at 37°C in an orbital shaker at 150 r.p.m. (New Brunswick, New Jersey) the bottles were centrifuged at 100 g for 4 min and the supernatant medium in each was replaced with 10 ml of growth medium. Each cell suspension was mixed thoroughly in a variable speed vibratory shaker (Fisons Limited), and a 35 µl sample was dispersed evenly in 5 ml growth medium in a 50 mm tissue culture dish (Nunclon Delta). The dishes were incubated for 7 days at 37°C in a humid atmosphere of 5% CO<sub>2</sub> in air. The cultures were then fixed with methanol and stained with Giemsa and the colonies in each dish were manually counted.

From the toxicity data 5 doses were selected for use in the transformation assays.

#### Cell Transformation Assay

Triplicate plates were used at each dose level of test compound and all experiments were repeated.

For all tests the procedure was described below.

Samples of BHK 21 C13 cells were harvested, suspended in MEMS with or without S-9 mix (5% v/v), and distributed among sterile plastic universal bottles as in the preliminary toxicity test.

The compound was dissolved in dimethylsulphoxide at the following concentrations:

200, 100, 50, 25 and 12.5 mg/ml

Triplicate cell suspensions, with and without S-9 mix, received 10  $\mu$ l samples of the test solutions.

The positive control compounds were N-methyl-N'-nitro-N-nitrosoguanidine (MNNG) for the S-9-free suspensions and 2-acetamidofluorene (2-AAF) for the S-9-supplemented suspensions. These were dissolved in dimethylsulphoxide in concentrations giving on dilution 1-100 in cell suspension the following levels:

MNNG:

0.3, 0.15, 0.075 and 0.375  $\mu$ g/ml

2-AAF:

600, 30υ, 150 and 75 μg/ml

Triplicate negative control suspensions received 10  $\mu l$  of dimethylsulphoxide.

Incubation of the suspensions, resuspension of the cells is growth medium and preparation of the cultures for estimation of cell survival were carried out as in the preliminary toxicity test.

After removal of the 35  $\mu$ l samples from each bottle for the measurement of toxicity, 625  $\mu$ l of 5% solution of agar (Difco Noble) in distilled water was mixed rapidly with the remaining suspension and the mixture poured into a 50 mm bacteriological Petri dish which was left at 4°C for 5 min to allow the agar to gel. The dishes were then incubated for 21 days at 37°C in sealed containers in a humid atmosphere of 5% CO<sub>2</sub> in air.

### Quantitation of Transformation

After a 7 day incubation period, the colonies growing in the tissue culture dishes were fixed for 30 min in methanol and stained with Giemsa. The colonies in each dish were counted manually and cell viability at each concentration of test compound was expressed as a percentage of the survivors in the negative control cultures, i.e. those dosed with dimethyl-sulphoxide only.

After 21 days incubation, the cells growing in soft agar were examined with a New Brunswick Biotran II Automatic Colony Counter (New Brunswick, New Jersey) and the number of transformants regarded as colonies with a diameter of > 0.22 mm (Bouck et al 1976), were counted. From the results of the transformation and survival assays the transformation frequency (number of transformed colonies/10<sup>5</sup> surviving cells) was calculated for each dose.

#### Evaluation

In the cell transformation test described by Styles (1977), a transformation frequency, at LC<sub>50</sub> of the compound under test, of 5 times that for the spontaneous transformation frequency is considered to indicate potential carcinogenicity. Using this method in tests on 120 compounds 91% of carcinogens and 97% of non-carcinogens were correctly distinguished, Purchase et al (1976). This is one criterion applied in the present study to indicate potential carcinogenicity of the test agent. In addition, if the test compound is relatively non-toxic or is insufficiently soluble to thieve a concentration resulting in 50% toxicity, a 2-fold increase in the absolute number of transformed colonies per dish at 2 doses is considered to indicate a positive response.

The clone of BHK 21 C13 cells used in these studies may in particular experiments give very low spontaneous transformation levels, i.e. 1 or 2. If this occurs the values obtained with the treated cells are related instead to the average spontaneous transformation rate (5 transformed colonies) found from a series of tests.

#### R'SULTS AND DISCUSSION

From the results of the toxicity test (Table 1), a top dose of 2000  $\mu g/ml$  was selected for the transformation assays. The addition of S-9 mix to the incubation mixture appeared to have only a slight effect on the toxicity of HE 1003, the mean LC<sub>50</sub> values in the presence and absence of S-9 mix being 1380 (Figures 1-3) and 1140  $\mu g/ml$  (Figures 5-7) respectively.

In the presence of S-9 mix, HE 1003, at the  $LC_{50}$ , caused 47.2, 13.2, 1.9 and 10.4 fold increases in transformation frequencies relative to that of the negative controls (Figures 1-4). These significant increases in the transformation frequencies at the  $LC_{50}$  were also accompanied by large increases in the absolute numbers of transformed colonies in cells treated with relatively non-toxic concentrations of HE 1003 (Tables 2-5).

In the absence of S-9 mix 3.1, 5.6 and 14 fold increases in transformation frequencies at the  $LC_{50}$  were measured (Figures 5-7). This, together with the increase in absolute numbers of transformed colonies in cells treated with HE 1003 below the  $LC_{50}$  level, indicates that the compound is a potential cell transforming agent (Tables 5-8).

The results obtained in this study indicate that HE 1003 should be regarded as a potential cell transforming agent.

# CONCLUSION

By the criterion used in this assay, HE 1003 showed evidence of cell transforming potential.

#### REFERENCES

- (1) Bouck, N and di Mayorca, G. (1976), Nature, 264, 722.
- (2) di Mayorca, G., Greenblatt, M., Trauthen, T., Soller, A. and Giordano, R. (1973), Proc. Nat. Acad. Sci., 70, 46.
- (3) Kao, F. and Puck, T.T. (1975), Genetics, 79, 343.
- (4) Purchase, I.F.H., Longstaff, E., Ashbey, J., Styles, J.A., Anderson, D., Lefeure, P.A. and Westwood, F.R. (1976), Nature, 264, 624.
- (5) Shin, S., Freedman, V.H., Risser, R. and Pollack, R. (1975), Proc. Nat. Acad. Sci., 72, 4435.
- (6) Styles, J. (1977), Brit. J. Cancer, 36, 558.

TABLE 1

# Cytotoxicity Test

Substance: Project No: 703863 Activation: None Contractor: Bayer Operator(s): Alan Poole Date of test: 19 August 1980

Liver preparation date: None Cell culture batch: 79.6 (4.1) Date survival assay counted: 26 August 1980

HE 1003

Survival as Survival Assay Colonies/plate Percentage of Negative Control Substance Quantity µg/ml Average 47 79 90, 68 1000 143.5 100 130, 157 97 162 165, 159 10 169.5 101 181, 158 1.0 103 171.5 179, 164 0.1 DMSO 100 167 170, 164 Control

#### TABL.

#### Cell Transfor

Project No: 703863

Substance: HE 1003

Contractor:

Bayer

Operator(s): Alan Poole, Nick Hunter, Catherine Green

Activation: Aroclor-induced Fischer Rat Liver preparation date: 13 August 1980

Date of test: 28 August 1980

Cell culture batch: 78.13 (2)

Date survival assay counted: 3 September 1980

Date transformation assay counted: 17 September 1980

Survival Assay Colonies/plate		Transformation Assay Colonies/plate		Viable Count	Transformed Colonies/105	Survival as Percentage of Negative
	Average		Average	x 10		Control
8, 0,	3.3	•, 4, 11	7.5	0.09	833	2
230,142,	161.3	67	67	4.61	145	85
212,206,	196.7	140,171,	153	5.62	272	104
159,195,	185.3	61, 97,	79	5.29	149	98
144,188,	173.7	63, 69, 78	70	4.96	141	92
0, 0,	0.3	0, 0.	1	0.0086	1163	0
13, 1,	5	0, 6,	4	0.14	286	3
8, 22, 37	22.3	2, 7,	5.7	0.64	89	12
216,170, 177	187.7	•, 3, 8	5.5	5.36	10	99
184,186,	189.3	5, 7,	4.7	5.41	9	100
	8, 0, 2 230,142, 112 212,206, 172 159,195, 202 144,188, 189  0, 0, 1 13, 1, 1 8, 22, 37 216,170, 177	Ryerage  8, 0, 3.3 2 230,142, 161.3 112 212,206, 196.7 172 159,195, 185.3 202 144,188, 173.7 189  0, 0, 0.3 1 13, 1, 5 1 8, 22, 22.3 37 216,170, 187.7 177	Survival Assay Colonies/plate  Average  8, 0, 3.3 e, 4, 2 230,142, 161.3 e, e, 112 212,206, 196.7 140,171, 172 148 159,195, 185.3 61, 97, 202 144,188, 173.7 63, 69, 189  0, 0, 0.3 0, 0. 1 3 13, 1, 5 0, 6, 1 8, 22, 22.3 2, 7, 37 8 216,170, 187.7 e, 3, 177  184,186, 189.3 5, 7,	Average	Survival Assay Colonies/plate    Average	Servival Assay   Colonies/plate   Count   Colonies/105

e Petri dish contaminated

#### Cell Transformation

Project No: 703863

HE 1003 Substance:

Contractor: Bayer

Operator(s): Alan Poole, Nick Hunter, Catherine Green

Activation: Aroclor-induced Fischer Rat Liver preparation date: 13 August 1980

Date of test: 8 September 1980

Cell culture batch: 78.13 (4)

Date survival assay counted: 19 September 1980

Date transformation assay counted: 29 September 1980

Substance Quantity µg/ml	Survival Assay Colonies/plate		Transformation Assay Colonies/plate		Viable Count	Transformed Colonies/10 <sup>5</sup> Viable Cells	Survival as Percentage or Negative
		Average		Average	× 10 <sup>4</sup>		Control
HE 1003							
2000	110,153	131.5	36,19,	27.5	3.76	73	34
1000	173,413,	299	19,39, 81	46.3	8.54	54	77
500	385,262,	307.7	25,38,	31.5	8.79	36	63
250	248,324,	309	22,34,	28.0	8.83	32	80
125	264,402, 293	319.7	6,59, 21	28.7	9.13	31	83
2-AAF							
600	2, 8,	4	1, 4,	3.0	0.11	273	1
300	3, 7,	3.7	4, 2,	3.0	0.11	273	1
150	8, 7, 9	8	8,10,	9.0	0.23	391	2
75	451,409, 272	377.3	12, 7,	9.5	10.78	9	98
DMSO Control	484,354, 321	386.3	5, 5,	5	11.04	5	100

e Petri dish contaminated

#### Cell Transformation

Project No: 703863 Contractor: Bayer

Operator(s): Alan Poole, Nick Hunter, Catherine Green

Date of test: 11 September 1980

Date survival assay counted: 22 September 1980
Date transformation assay counted: 2 October 1980

Substance: HE 1003

Activation: Aroclor-induced Fischer Rat Liver preparation date: 13 August 1980

Cell culture batch: 79.5 (12)

Substance Quantity	Survival Assay Colonies/plate		Transformation Assay Colonies/plate		Viable Count x 10 <sup>4</sup>	Transformed Colonies/105 Viable Cells	Survival as rercentage of Negative
µg/ml		Average		Average	X 10	The same	Control
HE 1003							N. Ma. all N
2000	0, 0,	0	0, 0,	0	0	0	0
1000	296,308, 278	294	20, •,	21.5	8.40	26	89
500	306,301	300	6, 4,	3.3	8.57	4	91
250	309,275,	284	6, 0,	3.0	8.11	4	86
125	306,298,	301.7	1, 6,	3.0	8,62	3	92
2-AAF							
600	48, 44,	48.0	7, 4,	4.3	1.37	31	15
300	104, 50,	95.0	0,12,	10.5	2.71	39	29
150	163,107,	133.0	41,40,	38.3	3.80	101	40
75	129 241,255, 256	250.7	0, 4,	1.3	7.16	2	76
DMSC Control	353,332, 301	328.7	4, 0,	7.5	9.39	8	100

e Petri dish contaminated

#### Cell Transformation

Project No: 703863 Contractor: Bayer

Contractor: Bayer
Operator(s): Alan Poole, Nick Hunter, Catherine Green

Date of test: 16 September 1980

Date survival assay counted: 26 September 1980
Date transformation assay counted: 7 October 1980

Substance: HE 1003

Activation: Aroclor-induced Fischer Rat Liver preparation date: 13 August 1980

Cell culture batch: 78.13 (3)

Subs_ance Quantity µg/ml	Survival Assay Colonies/plate		Transformation Assay Colonies/plate		Viable Count x 10 <sup>4</sup>	Transformed Colonies/105 Viable Cells	Survival as Percentage of Negative
		Average		Average	X 10		Control
HE 1003				V		the land	
2000	0, 0,	7.7	0, 0,	0	0.22	0	2
1000	229,291,	250.3	0,49,	49.0	7.15	69	67
500	208,223,	215.5	74,49, 49	57.3	6.16	93	58
250	219,326, 396	313.7	12, 0,	12.0	8.96	13	83
125	249,267, 317	277.7	25,31, 34	30.0	7.93	38	74
2-AAF							
600	1, 0,	0.5	1, 3,	3.7	0.01	3700	0
300	49, 1,	21.7	26, 4,	15.0	0.62	242	6
150	16, 46, 95	52.3	26,18,	22.0	1.50	147	14
75	253,372, 306	310.3	25,14, 12	17.0	8.87	19	83
DMSO Control	432,307,	374.3	5, 6,	5.0	10.70	5	100

e Petri dish contaminated

#### Cell Transformation

Project No: 703863

HE 1003 Substance:

Bayer Contractor:

Activation: None

Operator(s): Alan Poole, Nick Hunter, Catherine Green

Liver preparation date: None

Cell culture batch: 78.13 (4)

Date of test: 8 September 1980

Date survival assay counted: 22 September 1980

Date transformation assay counted: 29 September 1980

Substance Quantity µg/ml	Survival Assay Colonies/plate		Transformation Assay Colonies/plate		Viable Count x 10 <sup>4</sup>	Transformed Colonies/10 <sup>5</sup> Viable Cells	Survival as Percentage of Negative
		Average		Average	× 10		Control
HE 1003		VIII - I				4	
2000	1, 5,	3.7	1, 2,	1.0	0.11	91	1
1000	241,187,	205	8,12,	10.0	5.86	17	73
500	213,362,	283	28,25, 27	26.7	8.09	33	103
250	228,261,	241.3	14,17,	14.3	6.90	21	92
125	235,265,	243.3	6, 4,	5.0	6.95	7	93
MLNG							
0.3	0, 0,	0	0, 0,	0	0	0	0
0.15	65, 71,	45.3	18,17,	14.3	1.29	111	17
0.075	212,165,	200.3	14,10,	12.3	5.72	22	77
0.0375	242,277, 270	246.3	6, 4,	3.7	7.04	5	91
DMSO Control	316,236,	261.7	6,11, 15	10.7	7.48	14	100

e Petri dish contaminated

#### Call Transformation

Project No: 703863 Contractor: Bayer

r

Operator (s): Alan Poole, Nick Hunter, Catherine Green

Date of test: 11 September 1980

Date survival assay counted: 22 September 1980

Date transformation assay counted: 2 October 1980

Substance: HE 1003 Activation: None

Liver preparation date: None Cell culture batch: 79.5 (12)

Substance Co Quantity µg/ml	Survival Assay Colonies/plate		Transformation Assay Colonies/plate		Viable Count x 10 <sup>4</sup>	Transformed Colonies/105 Viable Cells	Survival as Percentage of Negative
		Average		Average	X 10		Control
HE 1003							
2000	0, 0,	0.7	0, •,	0	0.02	0	0
1000	180,228,	255.7	49,19,	34.0	7.3	47	83
500	359 354,248,	286.3	40,30,	36.0	8.18	44	93
250	257 294,245,	277.3	38 14,18,	17.7	7.92	22	90
125	293 417,267,	342.0	21 13, 4, 2	6.3	9.77	6	112
MNNG		1				168	15
0.3	35, 39, 63	45.7	16,26,	22	1.31	W.	
0.15	246,141,	206.7	11,16,	12.7	5.91	21	67
0.075	291,222,	253	3,17,	9.3	7.23	13	83
0.0375	246 313,282 261	285.3	•,11, 5	8.0	8.15	10	93
DMSO Control	354,291, 275	306.7	6, 4, 4	4.7	8.76	5	100

e Petri dish contaminated

TABLE 8

#### Cell Transformation

Substance: HE 1003

Liver preparation date: None

Cell culture batch: 78.13 (3)

Activation: None

Project No: 703863 Contractor: Bayer

or: Bayer (e): Alan Poole, Nick Hunter, Catherine Gree


Operator(s): Alan Poole, Nick Hunter, Catherine Green
Date of test: 16 September 1980

Date survival assay counted: 26 September 1980
Date transformation assay counted: 7 October 1980


Substance Quantity µg/ml	Surviva Colonie	Survival Assay Colonies/plate		Transformation Assay Colonies/plate		Transformed Colonies/10 <sup>5</sup> Viable Cells	Survival as Percentage of Negative
		Average		Average	× 10 <sup>4</sup>		Control
HE 1003							
2000	0, 0,	1.3	0, 0,	)	0.04	0	0
1000	0, 12, 73	28.3	38, •, 25	31.5	0.81	389	8
500	179, •,	179	23, •,	12.5	5.11	24	53
250	356,269,	282.3	78,62, 46	62.0	8.07	77	84
125	330,138, 150	206.0	64,39, 58	53.7	5.89	91	61
MNNG							
0.3	10, 38,	36.7	•, 3, 5	4.0	1.05	38	11
0.15	130,149,	148.7	12, 8,	15.7	4.25	37	44
0.075	252,208,	256.7	12,11,	15.7	7.33	21	76
0.0375	318,237, 277	277.3	12, 5,	8.3	7.92	11	82
DMSO Control	324,315, 375	338	3, 3,	3.3	9.66		100

e Petri dish contaminated


Survival and Transformation Frequency of BHK 21 C13 Cells Treated with HE 1003 in the Presence of S-9 Activation Graph of Results from Table 2

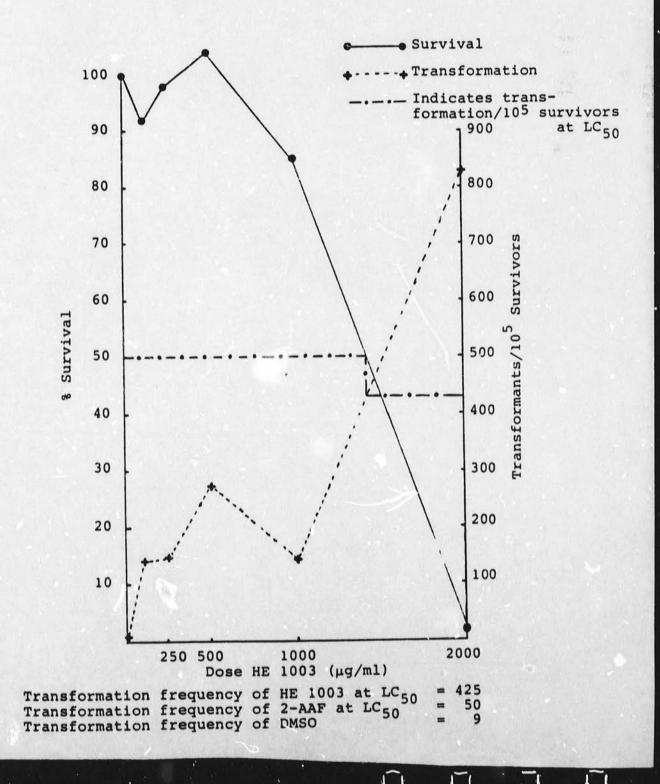


Survival and Transformation Frequency of BHK 21 C13 Cells Treated with HE 1003 in the Presence of S-9 Activation Graph of Results from Table 3

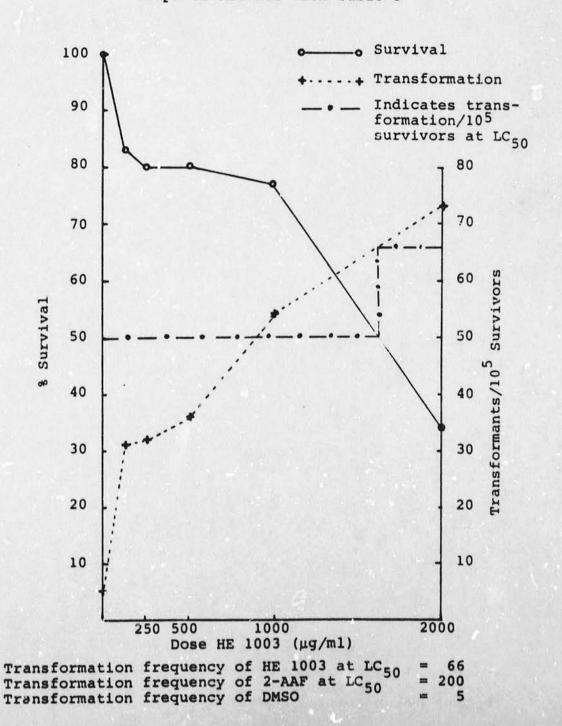


Survival and Transformation Frequency of BHK 21 Cl3 Cells Treated with HE 1003 in the Presence of S-9 Activation Graph of Results from Table 4

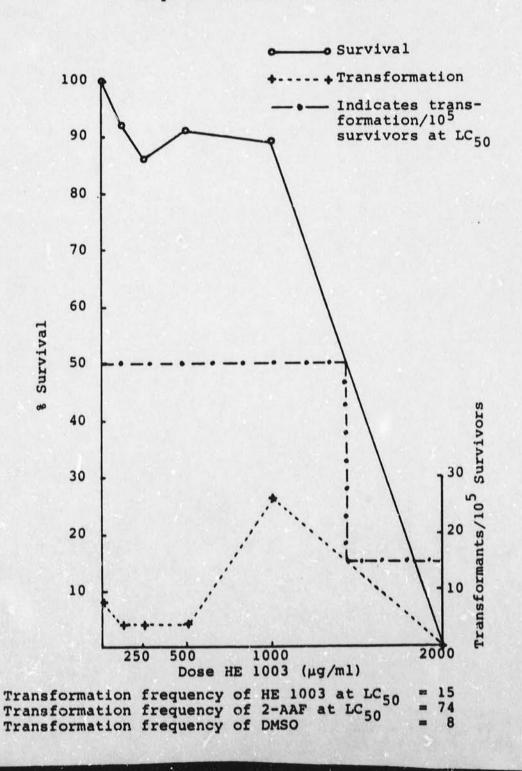


Survival and Transformation Frequency of BHK 21 C13 Cells Treated with HE 1003 in the Presence of S-9 Activation Graph of Results from Table 5

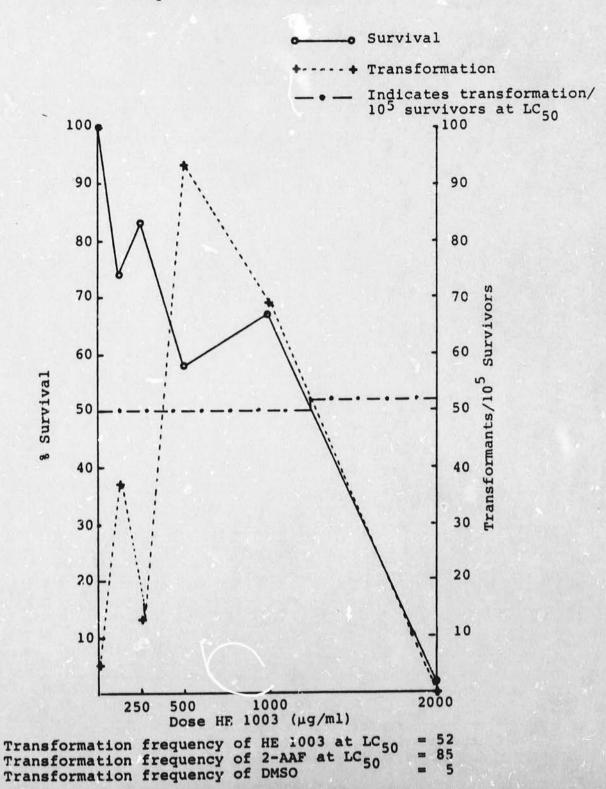


Survival and Transformation Frequency of BHK 21 C13 Cells Treated with HE 1003 in the Absence of S-9 Activation Graph of Results from Table 7

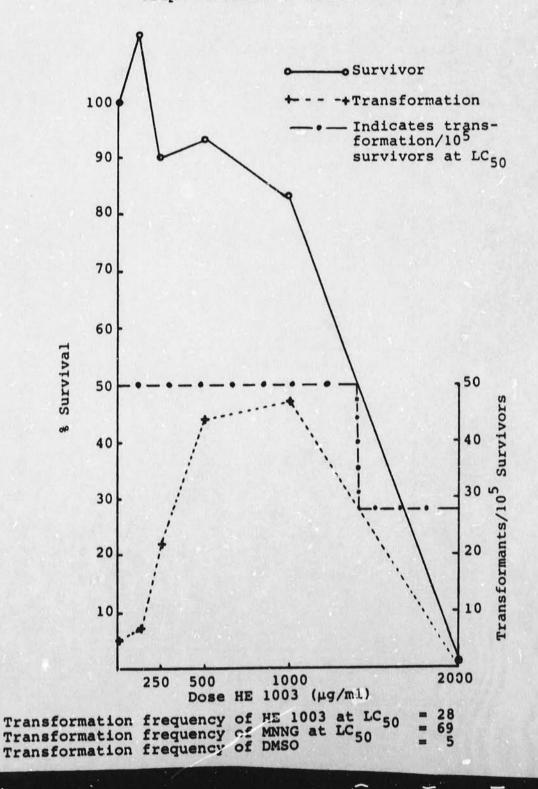

Survival and Transformation Frequency of BHK 21 C13 Cells Treated with HE 1003 in the Absence of S-9 Activation Graph of Results from Table 6

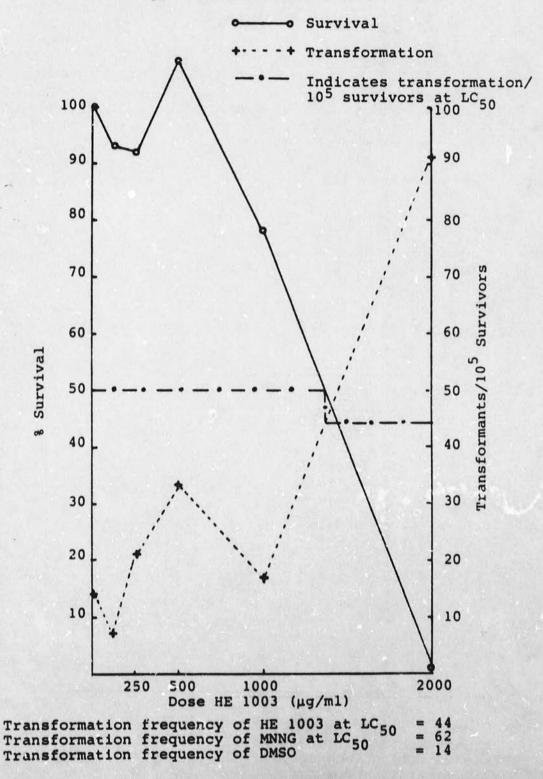

Survival and Transformation Frequency of BHK 21 C13 Cells Treated with HE 1003 in the Absence of S-9 Activation Graph of Results from Table 8

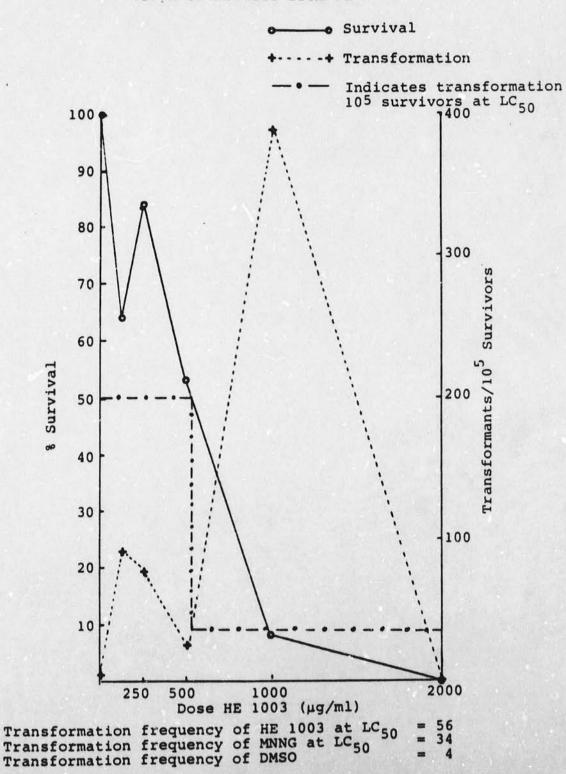

91 NOV 29 PM 3: 15

Charles of the care of the contract of the con

# CERTIFICATE OF AUTHENTICITY

THIS IS TO CERTIFY that the microimages appearing on this microfiche are accurate and complete reproductions of the records of U.S. Environmental Protection Agency documents as delivered in the regular course of business for microfilming.

Data produced 3 24 92 Marcia Lubalini
(Month) (Day) (Year) Camera Operator

Place Syracuse New York
(City) (State)

