

DOC COMPARISON STUDY
December 8, 1994
PAGE 1 OF 14

COMPARISON STUDY OF FIVE INSTRUMENTS MEASURING
DISSOLVED ORGANIC CARBON
FOR THE CHESAPEAKE BAY MONITORING PROGRAM

Betty A. Salley
Kevin Curling

December 8, 1994

ACKNOWLEDGEMENTS

This study would not have been feasible without the endeavors of the personnel at CBL, ODU, DHMH, and DCLS who analyze for dissolved organic carbon. They took time in their busy schedule to analyze the samples within the necessary time frame and forward the data.

TABLE OF CONTENTS

INTRODUCTION	4
HISTORY	4
STUDY DESIGN	5
INSTRUMENTS, METHODS, AND SAMPLE HANDLING	6
RESULTS AND DISCUSSION	8
CONCLUSIONS	12
REFERENCES	13
APPENDICES	
Appendix A:	Tables
Appendix B:	Figures
Appendix C:	Raw Data
Appendix D:	Instructions for analyzing samples
Appendix E:	Data and graphs from CSSP, Chesapeake Bay Program

**COMPARISON OF DISSOLVED ORGANIC CARBON INSTRUMENTS
AND EVALUATION OF THE EFFECT OF SALINITY**

INTRODUCTION

The Chesapeake Bay Monitoring Program (CBMP) was initiated by the U.S. Environmental Protection Agency, together with Maryland, Pennsylvania, Virginia and the District of Columbia in the summer of 1984. Water Quality monitoring has been performed on the Chesapeake Bay and its tributaries at least monthly since the program started. Dissolved Organic Carbon (DOC) has been measured from the beginning of the program.

The monitoring program was designed to develop a data base which would allow scientists 1) to determine trends in water quality over time and 2) to formulate a model of water quality processes. In order to provide continuity of data, a record of all methodology and instrument changes must be maintained.

The purpose of this study was two-fold:

- 1) Comparison of the DOC results from different instruments used in the program in order to provide a baseline for comparison of DOC data and historical data if these instruments are replaced.
- 2) Examination of salinity interference, if any, in of the results.

HISTORY

An interlaboratory study conducted at VIMS with two other laboratories suggested that DOC concentrations differed with the type of instrument used. Zimmermann (1991) suggested that these differences may have been due to salinity interference. In addition, samples which have been split between laboratories participating in the program have shown statistically significant differences. A study done in 1991 reported a significant difference in DOC when split samples were analyzed by two methods (Salley, 1991). The salinities of the samples were mesohaline to polyhaline.

In 1988, the international marine scientific community

became aware that a new DOC instrument was measuring much higher DOC concentrations than had been reported previously in open ocean water samples. This has led to intensive comparison of instruments and methodologies. An entire issue of Marine Chemistry (January, 1993) was devoted to the details of DOC measurement.

Due to this interest in and concern about DOC measurements, The Analytical Methods and Quality Assurance Workgroup (AMQAW) recommended that a comparison study be conducted between the DOC instruments involved in the Chesapeake Bay Monitoring Program at that time. The DOC study was also expanded to include an investigation to determine if the varying salinities of samples influenced or interfered with DOC measurements.

STUDY DESIGN

The study was designed to compare DOC recoveries by different instruments/methodologies. Since salinity possibly affected carbon recovery, the samples used for comparison were prepared in graduated salinities. Two carbon concentrations were used. The carbon concentrations chosen approximated the range encountered in the Chesapeake Bay. The five different instruments included in the study were used in the three mainstem Laboratories and two Tributary Laboratories. These instruments presently report DOC measurements to the Chesapeake Bay Monitoring Program (CBMP).

Sargasso Sea (SS) water was used for the high salinity diluent and Distilled and Deionized (DDI) water was used for the zero salinity diluent. The samples were prepared by combining the two waters to obtain the salinity required for the sample. For the low concentration DOC samples, the zero salinity carbon concentration was simply the residual carbon in the DDI water. Since the SS water contained more organic carbon than the DDI water, the carbon concentrations of the intermediate salinities were incrementally intermediate between those of the DDI and SS water. SS water and DDI water which had been separately adjusted to 8 mgC/L with glucose were used to prepare the high concentration DOC samples in the desired salinity range. Each sample was to have seven replicates analyzed. The replicates gave the statistical power for detection of differences in method and instrument variability.

INSTRUMENTS, METHODS, AND SAMPLE HANDLING

The procedures used in sample preparation, setting the blank, and standard curve were important parts of the results. Therefore, the method for each instrument is described in detail below.

1. Astro, Model 1850, Total Organic Carbon Analyzer

A 15 mL sample, pH <3, is purged with nitrogen and injected into the instrument. The sample is mixed with sodium persulfate and exposed to ultraviolet light. The resultant CO₂ is measured with a non-dispersive infrared cell. The standard curve is calculated using one standard and the instrument zero.

Low concentration samples were not analyzed on the Astro.

The study samples were acidified with 6 N H₂SO₄, transported, and stored at 4°C until analyzed.

2. Dohrmann, Model 180

A 5 mL sample is used. Phosphoric acid is added and the sample is sparged to remove the CO₂. Sodium persulfate is added and the sample is exposed to ultraviolet light. The instrument is zeroed with a millipore water blank. Five standards are analyzed and the data calculations are performed by the instrument.

The study samples were acidified with 6 N H₂SO₄, transported, and stored at 4°C until analyzed.

3. Oceanographic International Carbon Analyzer (OI) - Ampule method

A 5 ml sample, pH <3, is placed in an ampule and purged with ultrapure oxygen to remove the Dissolved Inorganic Carbon (DIC). One ml of saturated potassium persulfate and 200 µL of 10% phosphoric acid are added. The ampule is sealed and autoclaved at 130°C for four hours. The resultant CO₂ is carried through a Non-Dispersive Infrared Detector (NDIR) by nitrogen gas.

The NDIR is calibrated with blanks, standards and standard reference material (SRM) before the samples are analyzed. Check standards and spikes are interspersed throughout the analyses for internal quality control. The standards are reduced by linear regression and the intercept set to zero.

The study samples were acidified with 6 N H₂SO₄ and chilled to 4°C, for preservation and transport. On arrival, they were frozen until prepared for analysis. Since the pH was less than 3, no further acid was added before purging.

4. Oceanographic International Analytical, Model 700 TOC - Automated

The automated OI uses zero grade nitrogen as the carrier gas. All reagents (potassium persulfate and sodium phosphoric acid) are added automatically by the instrument. The CO₂ is purged from the digestion vessel with the nitrogen gas after phosphoric acid is added. In order to compensate for any possible salt interference in the oxidation, the volume of oxidant and the reaction time is increased when analyzing saline samples. Reagent blanks are run until instrument is stabilized. The standard curve is set with a one-point calibration where the standard concentration is 10 mgC/L. The zero point is set electronically. A scaling factor is calculated from this curve and used to calculate the sample values.

The study samples were not acidified, but frozen and transported. They remained frozen until analysis.

5. Shimadzu TOC 500, ASI-502, Automated

The Shimadzu method employs high temperature (680°C) combustion with a platinum catalyst. The carrier and sparge gas is zero-grade air. A sample, pH <3, is sparged for 6 minutes to remove DIC. An 80 µL sample is autoinjected into the TC port. The resultant carbon is oxidized to CO₂, dehumidified, and measured with a NDIR.

The instrument has an internal microprocessor. Each sample is injected at least three separate times. A coefficient of variation is calculated. If the coefficient is unacceptable, then instrument makes additional injections until the maximum of five injections is reached. An internal decision of which injections to use for the calculation is made by the microprocessor and the mean peak area, the standard deviation, and CV are printed out.

Instead of using the two point curve generated by the microprocessor, five internal standards are used with each set of analyses (18 samples) to calculate a linear regression with the intercept set at zero. Spikes, standards, and standard reference material are

interspersed throughout the analyses for quality control.

The study samples were acidified to pH <3 with HCl, stored at 4°C and analyzed within 30 days.

RESULTS AND DISCUSSION

Comparison of the instruments and any interference from salinity will be discussed separately. Since there was no means to determine the absolute true value for the carbon in the samples in this study, there was no single correct value. In order to provide a common baseline for comparing the instruments, a DDI blank sample was prepared and sent with the samples. The blank was the same water sample as the low carbon concentration, zero salinity sample, except that it was in a separate container. All concentrations used in the comparison have had this DDI water blank value subtracted from each sample value reported.

As previously noted, the Astro instrument did not analyze the low concentration samples. The laboratory personnel considered the instrument unable to correctly measure low levels of DOC in saline samples. Plans were made to replace both the Astro and Dohrmann instruments shortly after this study.

Comparison of Instruments and Data

The results from the four instruments that analyzed the low carbon concentration samples were very similar with the exception of one instrument. It was noted that the Dohrmann instrument analyzed each sample twice versus the seven replicates for the other instruments. In Figure 1, the Dohrmann values appear to differ from the others. When using the mean combined value for each sample from the Shimadzu and OI instruments, the maximum and minimum differences from this mean by the Dohrmann is 1.32 and 0.19 mgC/L respectively. The average difference between the Dohrmann and the mean of the other instruments' values is 0.68 mgC/L.

When the data from the Shimadzu and the two OI instruments were compared, there was a maximum difference between the nine low carbon concentration samples of 0.35 mgC/L and a minimum difference of 0.02 mgC/L. In general, the difference increased slightly as the salinity increased. In most instances the values of the three instruments were within one standard deviation of each other. The minimal concentration of these samples does not encourage great

variability. If the blank sample value had not been subtracted from each sample mean, the difference would have been greater.

The high carbon concentration samples analyzed by the five instruments showed wider variability. Because the concentration of samples was 8 mgC/L or higher, there was a greater potential for differences in the absolute recovery than there was when comparing the absolute recoveries of the low carbon concentration samples. All five instruments' high carbon concentration determinations are plotted in Figure 2. The Dohrmann and Astro sample values showed a greater variability than the Shimadzu and the two OI results. In order to clarify the graph, two additional plots were made: Figure 3 which displays only the high carbon concentration determinations from the Shimadzu, OI ampule and OI automated instruments and Figure 4 which displays the mean of these three determinations plotted with the Dohrmann and Astro results. Assuming that the sample values ranged from 8 to 9 mgC/L, increasing with salinity, these two instruments were generally within 1 mgC/L of the mean value of the other instruments. As with the low carbon samples, the Dohrmann only analyzed two replicates as opposed to seven or more by the other instruments.

The standard deviation of the replicates indicates the variability within each instrument. Precision is not an indicator of accuracy, but it is indicative of the instrument's quality control. Using only the high concentration samples to illustrate, all values are in mgC/L;

- 1.) The Astro standard deviation varied from 0.1 to 0.59.
- 2.) The Dohrmann with two replicates can not be compared.
- 3.) The OI ampule standard deviation varied from 0.06 to 0.17.
- 4.) The OI automated standard deviation varied from 0.03 to 0.14.
- 5.) The Shimadzu standard deviation varied from 0.05 to 0.18.

When the results from the Shimadzu, OI ampule, and OI automated instrument were compared, the high carbon samples of varying salinities showed as little variability as the low carbon samples. The maximum difference was 0.67 mgC/L and the minimum difference was 0.12 mgC/L. The average difference for the three instruments was 0.32 mgC/L. This is within the Upper Control Limit for precision for some DOC instruments.

Influence of Salinity

Since the Shimadzu, OI ampule, and OI automated data were so comparable, their data were used to access the influence of salinity on the DOC results. The samples had been prepared to allow close examination of the salinity range from 5 to 14 ppt. Previous studies indicated that salinities in this range influenced DOC results. No salinity influence would give a straight line regression of the sample values and salinity for each set. On examination of Figures 1 and 3, each instrument showed some slight variability in the area of 5 to 14 ppt salinity. In order to superimpose the data from the three instruments, the sample set for low carbon concentration was set to zero for the zero salinity (Figure 5) and 8.0 mgC/L for the first high level sample (Figure 6). The rest of the sample means for each set were adjusted accordingly.

This does not address the possibility of a continual salinity effect, which would give a linear regression, but it allows comparison of the methodologies in analyzing saline samples.

Standards, Blanks, and Curves

Although this study did not include laboratory calibration of instruments and data reduction, the final results were very dependent on these factors. In addition, the problem of carbon contamination was of great interest and was very much a problem in analyzing low level carbon samples. Known sources of carbon contamination include:

- C1.) The acid added to lower the pH of samples and standards to <3 before sparging off the inorganic carbon. This is referred to as the Acid Blank.
- C2.) The water to prepare the standards. This is referred to as the Standard Blank.
- C3.) Carbon in the instrument through which the sample passes. This is referred to as the Instrument Blank.

The samples and standards during analysis are equally contaminated by carbon from sources C1. and C3. Only the standards contain carbon contamination from the Standard Blank (C2.).

In addition to contamination, methods of data reduction for the samples need to be evaluated. Examples of some methods are given below:

- S1.) Some instruments use an electronic zero, but not a standard blank. Using this electronic zero and a single high standard, a two point curve is regressed. The sample values are calculated from this curve.

S2.) In some instruments a set of standards containing from four to ten separate values is analyzed and a linear regression calculated using the results. The sample concentrations values are calculated using this curve.

a. This regression is allowed to set its own intercept which is generally above zero and no further adjustment is made. This assumes that any carbon found in the standards is also found in the samples.

b. Dividing the sample peak area by the slope of the regression could be used to calculate the sample values. This sets the curve to zero and assumes that any carbon found in the standards is not found in the samples.

When the method described in S1. is used, the electronic zero assumes no carbon in the Acid or Standard Blank; however, it makes allowance for the Instrument Blank. Therefore, the resultant standard curve would be slightly skewed. The values close to zero may be higher than the true value and the values in the upper range may be slightly lower than the true value due to the presence of carbon in the standard water for which no correction can be made.

When the method described in S2a. is used, there is no allowance in the calculation for the carbon in the Acid Blank, Standard Blank, and Instrument Blank; thus, these would be included in the resultant standard curve. Consequently, a lower sample result than the true value would be obtained.

When the S2b. method is used, the Acid Blank and Instrument Blank are fully considered, but the carbon in the standard water would be included in the resultant standard curve, thus giving sample concentrations which are higher than the true value.

In addition, other factors exist which can result in a bias in very low level carbon analysis. For instance, the standard diluent typically has some carbon contamination present which is difficult to remove and should be considered in the calculations. As the level of carbon in the sample increases, this contamination assumes less importance. Most DDI water contains no more than 0.30 mgC/L, when only polished with Deionization cartridges. When any water (samples, reagents, standards, etc...) is exposed to the atmosphere, it collects carbon. Further consideration is given to these biases in Sharp et al. (1994).

CONCLUSIONS

Comparison of the dissolved organic carbon results from the five methodologies/instruments used for the Chesapeake Monitoring Program demonstrated that the consensus was good between all instruments and excellent between three instruments. There was no need to develop a correction factor between instruments. Any problems with the two instruments which yielded the greatest variability was probably due to the age of the instruments and their detectors. That instrument variability has been eliminated as those two instruments have been replaced by high temperature combustion instruments which performed well in this comparison.

Salinity may present a bias in wet oxidation methods, but the instruments in this study were modified in their reagent concentration and reaction time to compensate for salinity. The motion linearity of the curves with the adjusted means for examination of salinity influence is probably an artifact of the addition measurements in that range. It could be argued from these curves that salinity does affect DOC measurements, but if so, the influence is so slight that it can be ignored.

Attention to calibration, blanks, and standards is more important than difference in recovery of dissolved organic carbon by the instruments. All data reported must have an explanation of how this was handled for future data users. However, for the concentrations of carbon found in the Chesapeake Bay Monitoring Program, errors induced in measurement by data reduction are probably less than those normally associated with intralaboratory variability since all the laboratories involved share a common method of data generation.

REFERENCES

- Benner, R. and Strom, M., 1992. A critical evaluation of the analytical blank associated with DOC measurement by high-temperature catalytic oxidation. *Marine Chemistry*, 41:153-160.
- Hedges, J. I., Bergamaschi, B. A. and Benner, R. 1993. Comparative analyses of DOC and DON in natural water. *Marine Chemistry*, 41: 121-134.
- Perdue, E.M. and Mantoura, F. 1993. Mechanism Subgroup Report. *Marine Chemistry*, 41:51-60.
- Salley, B. A., Curling, K. & Neilson, B, 1991. A comparison of two DOC methods. Report to Analytical Methods Quality Assurance Workgroup, Chesapeake Bay Program.
- Sharp, J.H., 1973. Total organic carbon in seawater - comparison of measurements using persulfate oxidation and high temperature combustion. *Marine Chemistry*, 1:211-229.
- Sharp, J.H., Suzuki, Y. & Munday, W.L. 1988. Intercomparison of dissolved organic carbon analyses in estuarine and coastal waters of the North Atlantic Ocean. *EOS*, 69: 1134.
- Sharp, J.H., Benner, R., Bennett, L., Carlson, C.A., Fitzwater, S.E., Peltzer, E.T., & Tupas, L.M. 1994. Analyses of Dissolved Organic Carbon in Seawater: The JGOFS EQPAC Methods Comparison, in press.
- Sharp, J.H. 1993. The Dissolved Organic Carbon Controversy: An Update. *Oceanography*, 6(2): 45-50.
- Shimadzu Corporation, 1989. Instruction Manual, Total Organic carbon Analyzer, Auto Sample Injection, Model ASI-502, CM 393-070 and 638,90887. Kyoto, Japan.
- Suzuki, Y. 1993. On the Measurement of DOC and DON in seawater. *Marine Chemistry*, 41: 287.
- Tupas, L.M., Popp, B.N., & Karl, D.M. 1994. Dissolved Organic Carbon in Oligotrophic Waters: Experiments on Sample Preservation, Storage, and Analysis. *Marine Chemistry*, 45: 207-216.
- US Environmental Protection Agency, 1983. Methods for Chemical Analysis of Water and Wastes. Environmental Research Center, Cincinnati, Ohio.

US Environmental Protection Agency, 1990-1991. Chesapeake Bay Coordinated Split Sample Program Annual Report. Chesapeake Bay Program Office, Annapolis, Maryland.

Wangersky, P. J., 1993. Dissolved organic carbon methods: a critical review. Marine Chemistry, 41: 61-74.

Williams, P. M. and E. R. M. Druffel, 1988. Dissolved Organic Matter in the Ocean: Comments on a Controversy. Oceanography, 1:1, 14-17

Williams, P., Bauer, J., 1993. DOC subgroup report. Marine Chemistry, 41:11-21.

Zimmermann, C. F., 1991. Estuarine Nutrient Analyses: a comparison of sample handling techniques and analysis of carbon, nitrogen, phosphorus and Chlorophyll a. U.S.EPA Environmental Monitoring Systems Laboratory, Cincinnati, OH 45268. 12 pages.

APPENDICES

- A. Tables
- B. Figures
- C. Raw Data
- D. Instruction for Analyses
- E. Coordinated Split Sample Data from Chesapeake Bay Program

APPENDIX A

Tables

Contents:

Table A1. Instrument Methodology

Table A2. Means of DOC data by Salinity for Low Concentration

Table A3. Means of DOC data by Salinity for High
Concentration

Table A4. Adjusted Means of DOC data for Low Concentration

Table A5. Adjusted Means of DOC data for Low Concentration

TABLE 1

INSTRUMENT METHODOLOGY

ITEM	ASTRO	DOHRMANN	OI AMPULE	OI AUTO	SHIMADZU
SAMPLE SIZE	15mL	5 mL	5 mL	1.103mL	80µL
ACID ADDED	H2SO4	H2SO4	H2SO4	H3PO4 added internally	H2SO4
METHOD	Sodium Persulfate/ ultraviolet/ NDIR detector	SAME AS ASTRO	Sodium Persulfate/ 10% H3PO4 autoclaved	Sodium persulfate/ H3PO4 100°C	Platinum Catalyst 680°C
CALIBRATION	ONE POINT CURVE AND INSTRUMENT BLANK	TWO POINT CURVE AND MILLIPORE WATER FOR BLANK	ELEVEN POINTS PLUS MILLIPORE BLANK	ONE POINT CALIBRATION WITH INSTRUMENT BLANK	FIVE POINT CALIBRATION
GAS	Nitrogen purge	Nitrogen purge	Nitrogen purge	Zero Grade N ₂	Zero Grade Air
DATA REDUCTION	Regressed on calibration curve	Calculation by Instrument	Double Regression with intercept at zero	Instrument calculates a scaling factor for samples mVs	Slope of standards calculated & area count of samples divided by slope

TABLE 2.

MEANS OF DOC DATA IN MG/L BY SALINITY FOR LOW CONCENTRATIONS
corrected for blank value

SAMPLE ID	SALINITY	DORHMANN	OI AMPULE	OI AUTOMATED	SHIMADZU TOC 500
L7	0.03	0.90	0.02	-0.04	0.06
L1	5.20	1.54	0.21	0.19	0.26
L2	7.18	0.92	0.19	0.25	0.31
L9	9.23	0.58	0.19	0.32	0.27
L5	11.36	0.87	0.32	0.33	0.34
L6	12.99	1.00	0.31	0.35	0.44
L4	14.97	0.65	0.44	0.44	0.50
L3	24.53	1.15	0.77	0.68	0.84
L8	34.66	2.14	1.00	0.88	1.23

TABLE 3.

MEANS OF DOC DATA IN MG/L BY SALINITY FOR HIGH CONCENTRATIONS
corrected for blank value

SAMPLE ID	SALINITY	ASTRO	DORHMANN	OI AMPULE	OI AUTO	SHIMADZU TOC 500
H8	0.06	8.78	8.72	8.12	7.80	8.03
H6	5.15	9.65	9.36	7.91	7.99	8.14
H3	7.25	9.39	7.37	8.53	8.14	8.27
H2	9.22	9.31	9.38	8.28	8.23	8.16
H5	11.27	9.90	9.33	8.11	8.25	8.41
H1	13.40	8.54	8.29	8.29	8.31	8.43
H4	15.40	7.97	9.65	8.53	8.35	8.37
H7	23.72	8.23	9.36	8.76	8.42	8.83
H9	35.50	9.39	9.96	9.35	8.68	9.19

TABLE 4.
 ADJUSTED MEANS OF DOC DATA FOR LOW CONCENTRATIONS
 BY SUBTRACTION OF ZERO SALINITY VALUE
 IN MG C/L

SAMPLE ID	SALINITY	OI AMPULE	OI AUTOMATED	SHIMADZU TOC 500
L7	0.03	0.00	0.00	0.00
L1	5.20	0.19	0.23	0.20
L2	7.18	0.17	0.29	0.25
L9	9.23	0.17	0.36	0.21
L5	11.36	0.30	0.37	0.28
L6	12.99	0.29	0.39	0.38
L4	14.97	0.42	0.44	0.44
L3	24.53	0.75	0.72	0.78
L8	34.66	0.98	0.92	1.17

TABLE 5.
 ADJUSTED MEANS OF DOC DATA FOR HIGH CONCENTRATIONS
 TO 8.0 MG C/L FOR ZERO SALINITY VALUE

SAMPLE ID	SALINITY	OI AMPULE	OI AUTO	SHIMADZU TOC 500
H8	0.06	8.00	8.00	8.00
H6	5.15	7.79	8.19	8.11
H3	7.25	8.41	8.34	8.24
H2	9.22	8.16	8.43	8.13
H5	11.27	7.99	8.45	8.38
H1	13.40	8.17	8.51	8.40
H4	15.40	8.41	8.55	8.34
H7	23.72	8.64	8.62	8.80
H9	35.50	9.23	8.88	9.16

Appendix B

Figures

Figure 1.

Figure 2.

Figure 3.

Figure 4.

Figure 5.

Figure 6.

DOC Comparison

LOW CONCENTRATION


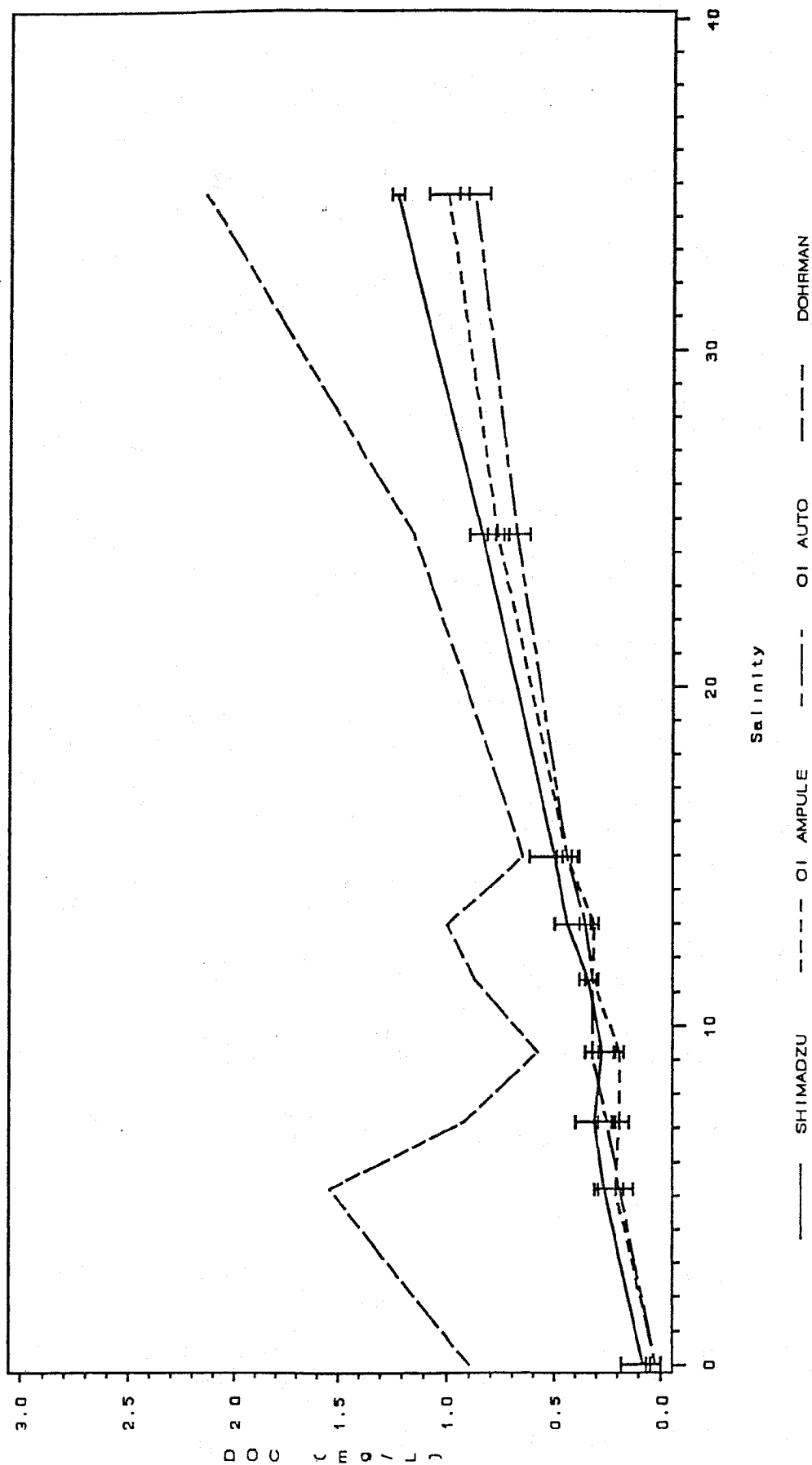


Figure 1.

DOC Comparison

HIGH CONCENTRATION

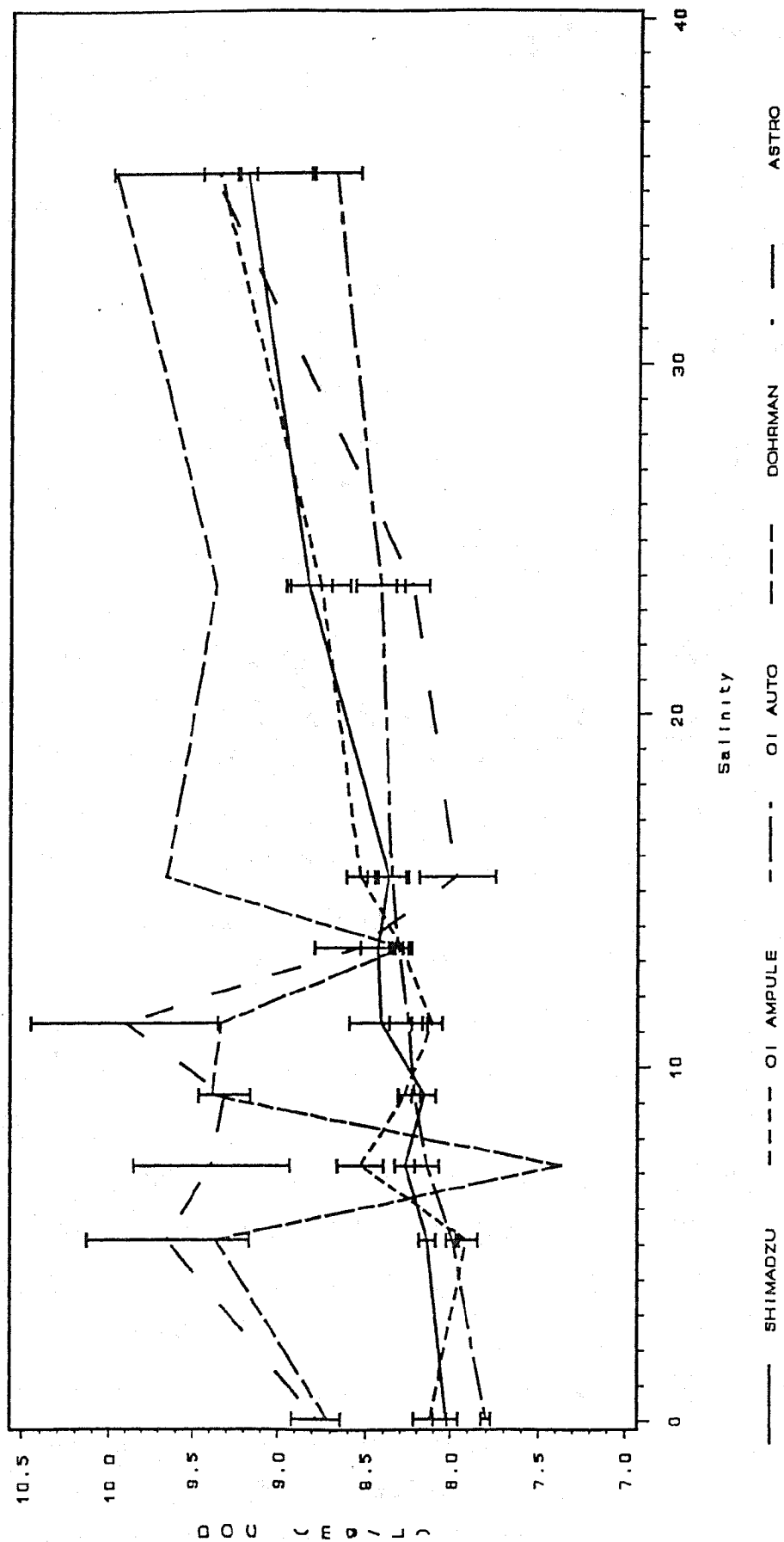

Mean of SHIMADZU, OI AUTO, OI AMPULE

Figure 2.

DOC Comparison

HIGH CONCENTRATION


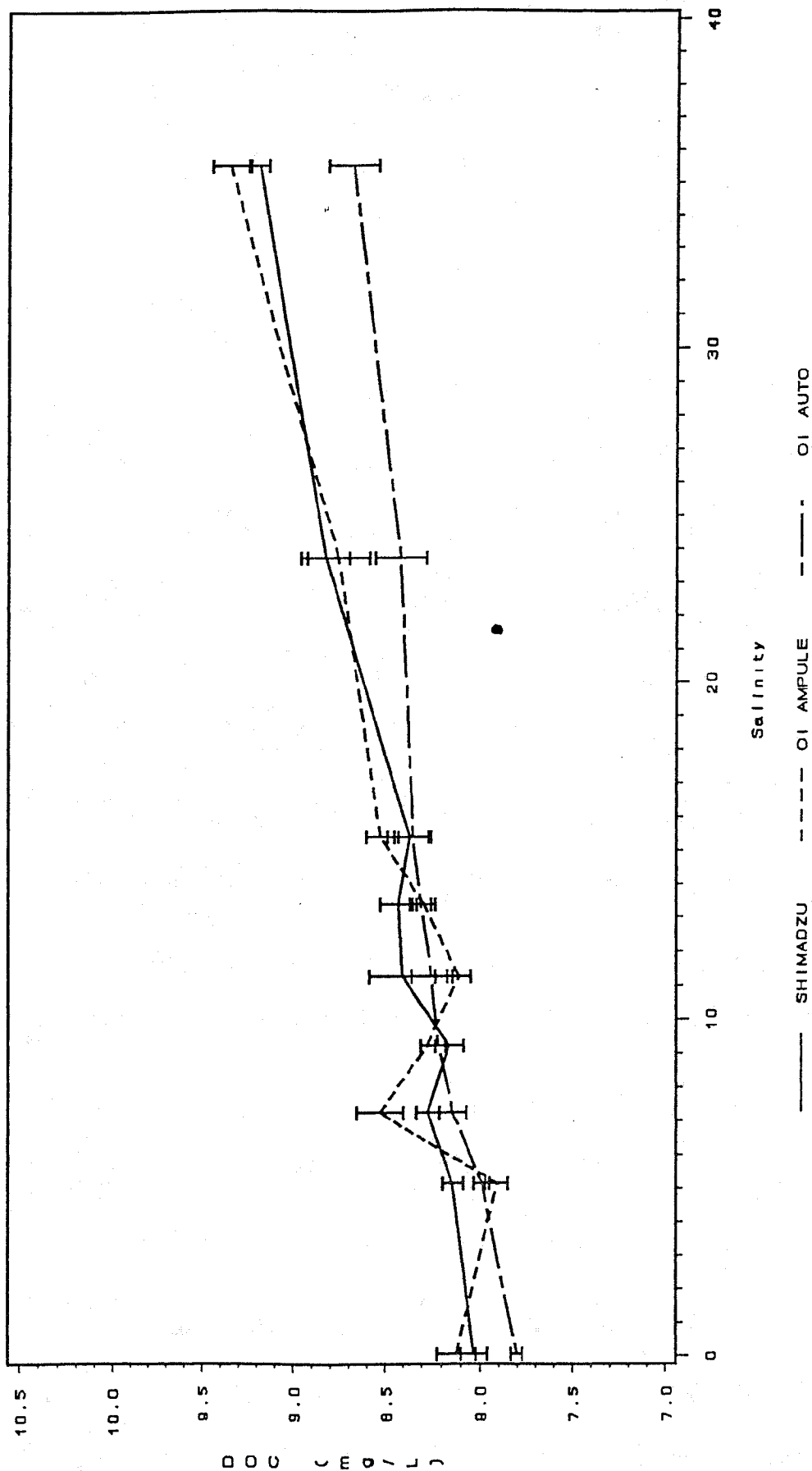


Figure 3.

DOC Comparison

HIGH CONCENTRATION

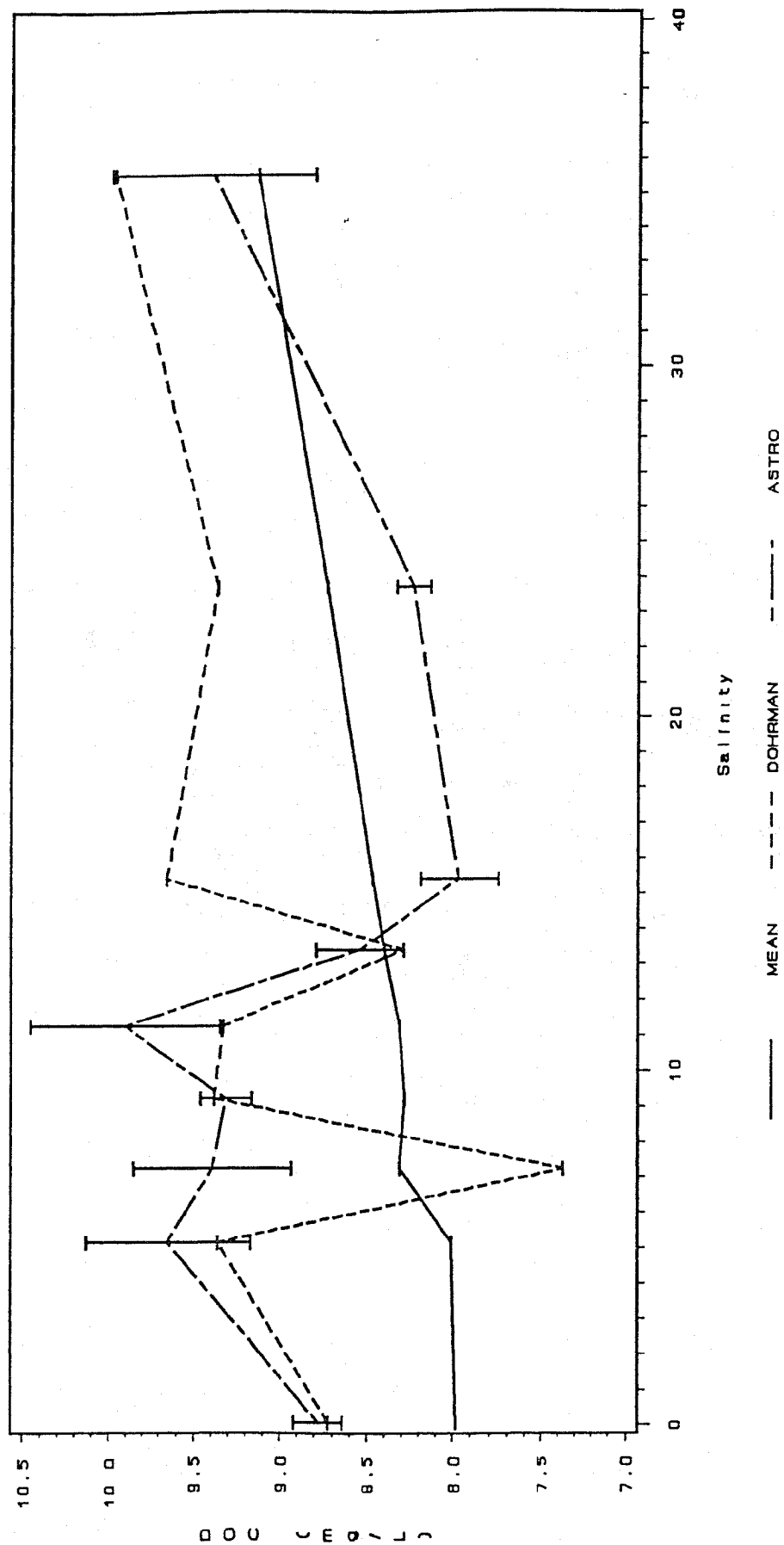

Mean of SHIMADZU, OI AUTO, OI AMPULE

Figure 4.

DOC Comparison

ADJUSTED LOW CONCENTRATION

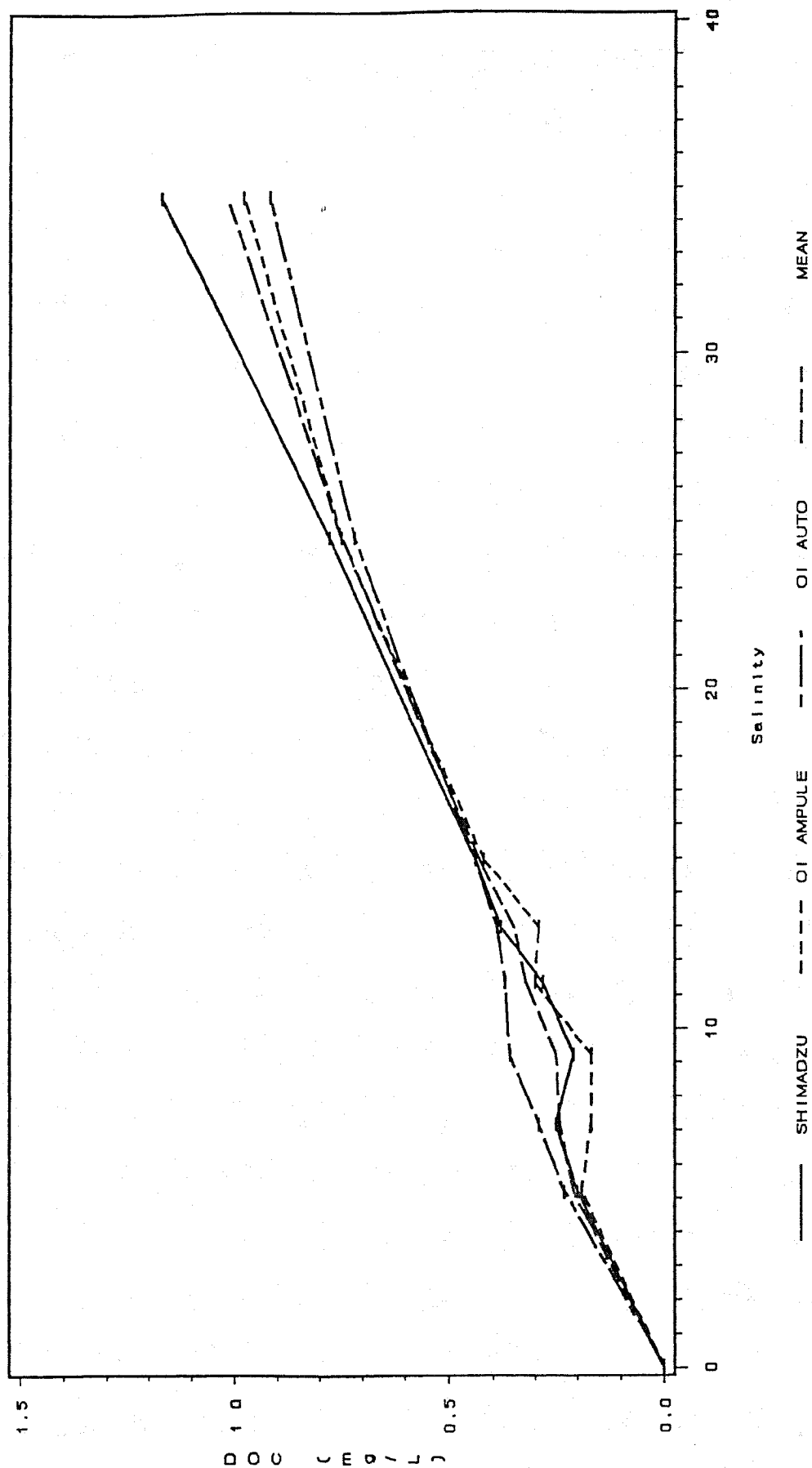


Figure 5

DOC Comparison

ADJUSTED HIGH CONCENTRATION

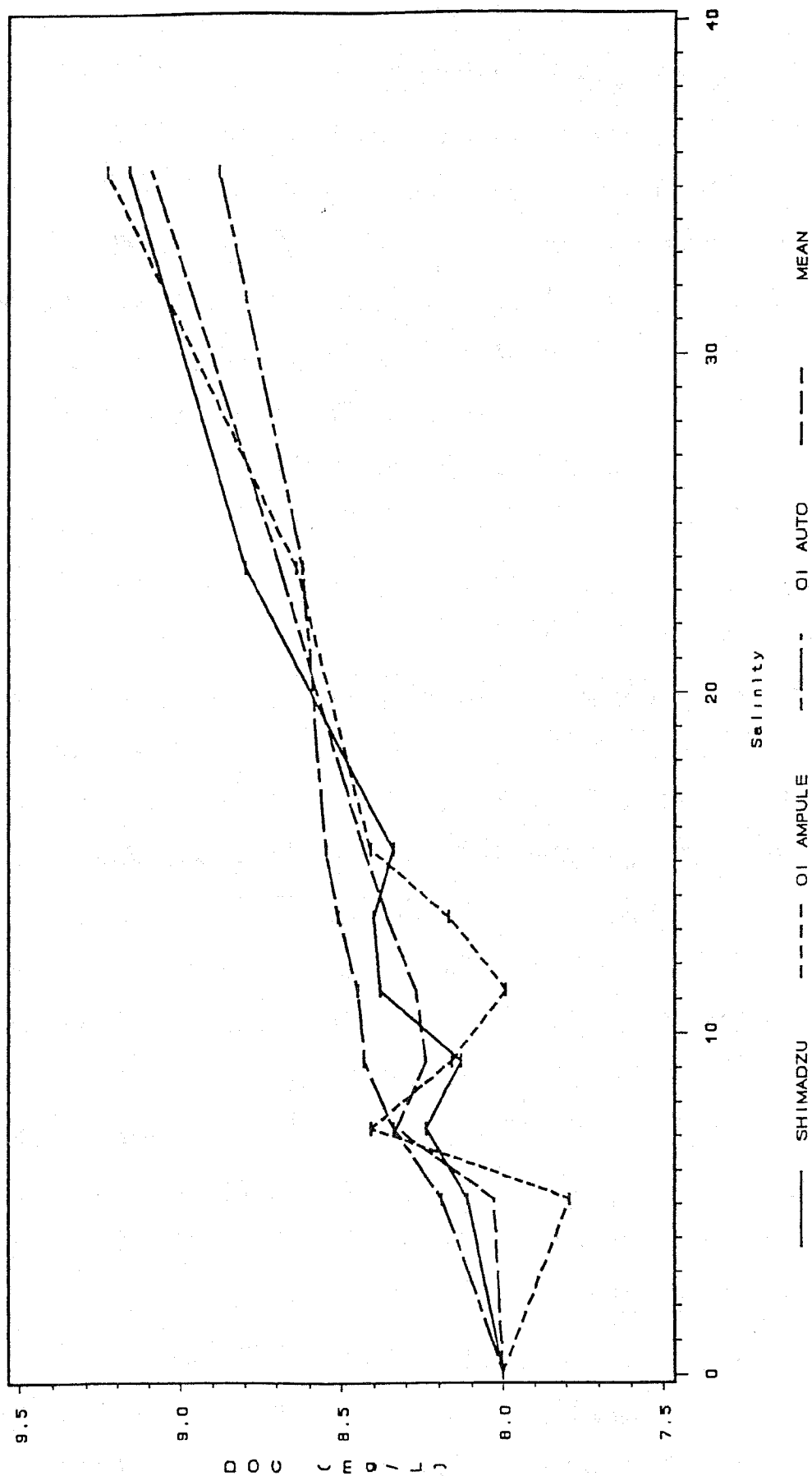


Figure 6.

APPENDIX C
Raw Data

1. Astro
2. Dohrmann
3. Oceanographic International, Ampule
4. Oceanographic International, Automated
5. Shimadzu

DHMH - Astro

VIMS DOC COMPARISON STUDY

in mg/L

DI H₂O BLANK

REP 1 0.40
2 0.42
3 0.40
4 0.41
5 0.42
6 0.45
7 0.41
MEAN 0.42
STD .016

	H1	H2	H3	H4	H5	H6	H7	H8	H9
REP 1	8.76	9.49	9.54	8.36	9.42	11.03	9.87	9.15	9.28
2	9.27	9.90	9.43	8.76	10.56	9.89	8.54	9.23	9.93
3	9.19	9.91	10.36	8.48	10.51	10.31	8.76	9.25	9.24
4	9.13	9.69	9.88	8.34	9.87	10.09	8.63	9.15	9.17
5	8.75	9.63	9.34	8.10	10.73	9.56	8.60	9.02	9.37
6	8.57	9.76	9.64	8.31	10.81	9.75	8.79	9.47	10.68
7	9.07		10.50		9.87			9.13	9.28
MEAN	8.96	9.73	9.81	8.39	10.32	10.07	8.65	9.20	9.61
STD	0.247	0.148	0.457	0.218	0.549	0.483	0.099	0.140	0.588
minus blank	8.54	9.31	9.39	7.97	9.90	9.65	8.23	8.78	9.39

SPIKES

	ACTUAL	ORIGINAL	EXPECTED
H4	13.76	9.81	14.81
H2	14.95	9.73	14.81
H5	14.57	10.32	15.32
H7	13.47	8.65	13.65

SRM True Value = 8.2 mg/L

1 8.62
2 8.79
3 8.63

DCLS - Dorchmann

VIMS DOC COMPARISON STUDY
in mg/L

	L1	L2	L3	L4	L5	L6	L7	L8	L9
REP 1	1.633	0.863	1.334	0.383	0.949	1.046	0.636	2.747	0.755
2	1.448	0.984	0.966	0.925	0.798	0.946	1.155	1.536	0.402
MEAN	1.541	0.924	1.150	0.654	0.874	0.996	0.896	2.142	0.579
minus									
blank	1.437	0.820	1.046	0.550	0.770	0.892	0.792	2.038	0.475

	H1	H2	H3	H4	H5	H6	H7	H8	H9
REP1	7.496	9.448	8.058	9.412	9.781	9.394	9.250	8.574	9.870
2	9.299	9.314	6.886	10.100	9.083	9.531	9.452	9.074	10.260
MEAN	8.398	9.381	7.472	9.756	9.432	9.463	9.463	8.824	10.065
minus									
blank	8.294	9.377	7.368	9.652	9.328	9.359	9.359	8.720	9.961

Blank	
Rep 1	0.113
Rep 2	0.095
Mean	0.104

ODU - Oceanographic International, Ampule

VIMS DOC COMPARISON STUDY
 in mg/L

	L1	L2	L3	L4	L5	L6	L7	L8	L9
REP 1	0.58	0.42	1.06	0.63	0.53	0.52	0.24	1.36	0.36
2	0.44	0.35	0.96	0.61	0.51	0.49	0.21	1.12	0.34
3	0.35	0.33	0.92	0.62	0.51	0.47	0.19	1.13	0.39
4	0.34	0.40	0.92	0.59	0.50	0.47	0.18	1.15	0.35
5	0.35	0.35	0.93	0.59	0.47	0.46	0.18	1.14	0.35
6	0.34	0.35	0.91	0.61	0.49	0.47	0.17	1.18	0.35
7	0.35	0.30	0.91	0.63	0.45	0.47	0.17	1.09	0.37
DUP 1	0.36								
2	0.33								
3	0.40								
MEAN	0.38	0.36	0.94	0.61	0.49	0.48	0.19	1.17	0.36
STD	0.077	0.041	0.054	0.017	0.027	0.020	0.025	0.089	0.017
minus blank	0.21	0.19	0.77	0.44	0.32	0.31	0.02	1.00	0.19

DI H2O

BLANK

REP 1	0.18
2	0.15
3	0.15
4	0.15
5	0.19
6	0.19
7	0.18

MEAN	0.17
STD	.019

ODU(2)

	H1	H2	H3	H4	H5	H6	H7	H8	H9
REP 1	8.46	8.39	8.42	8.57	8.21	8.07	9.27	8.18	9.55
2	8.38	8.45	8.64	8.66	8.22	8.07	9.19	8.31	9.56
3	8.37	8.55	8.67	8.67	8.25	8.01	8.85	8.26	9.44
4	8.54	8.46	8.81	8.74	8.27	8.05	8.86	8.31	9.31
5	8.51	8.43	8.74	8.79	8.29	8.08	8.86	8.48	9.58
6	8.49	8.42	8.76	8.71	8.31	8.19	8.74	8.21	9.56
7	8.48	8.43	8.77	8.78	8.40	8.10	8.99	8.28	9.59
DUP 1			8.82				8.84		9.58
2							8.82		
3							8.89		
MEAN	8.46	8.45	8.70	8.70	8.28	8.08	8.93	8.29	9.52
STD	.064	.051	.131	.077	.064	.056	.170	.097	.098
minus blank	8.29	8.28	8.53	8.53	8.11	7.91	8.76	8.12	9.35

SPIKES

	ACTUAL	ORIGINAL	EXPECTED
L2	4.49	0.36	4.36
L8	5.14	1.17	5.17
H4	12.53	8.70	12.70
H6	12.68	8.08	12.08

	EPA	
	3.06	9.18
		6.14
1	2.94	9.33
2	2.83	9.36
3	2.87	9.37
4	2.86	9.41

MEAN	2.88	9.37	6.13
STD	3.19	9.49	6.27

CBL - Oceanographic International, Automated

VIMS DOC COMPARISON STUDY
in mg/L

	L1	L2	L3	L4	L5	L6	L7	L8	L9
REP 1	0.39	0.51	0.88	0.71	0.55	0.50	0.13	1.03	0.52
2	0.34	0.39	0.98	0.70	0.50	0.51	0.13	1.17	0.48
3	0.39	0.40	0.79	0.60	0.47	0.53	0.16	1.10	0.54
4	0.36	0.42	0.80	0.58	0.49	0.58	0.12	1.02	0.51
5	0.41	0.38	0.88	0.60	0.49	0.58	0.15	1.09	0.52
6	0.34	0.46	0.88	0.58	0.49	0.52	0.20	1.05	0.46
7	0.37	0.43	0.82	0.57	0.51	0.52	0.15	0.90	0.48
DUP 1		0.43			0.54		0.11	1.07	
2								1.12	
3								1.04	

MEAN	0.37	0.43	0.86	0.62	0.51	0.53	0.14	1.06	0.50
STD	0.02	0.04	0.06	0.05	0.03	0.03	0.03	0.07	0.03
minus									
blank	0.19	0.25	0.68	0.44	0.33	0.35	-0.04	0.88	0.32

DDI H2O BLANK

REP 1	0.14
2	0.12
3	0.18
4	0.12
5	0.11
6	0.41
7	0.21

DUP 1	0.16
2	0.08
3	0.07
MEAN	0.18

CBL (2)

	H1	H2	H3	H4	H5	H6	H7	H8	H9
REP 1	8.50	8.35	8.31	8.41	8.34	8.12	8.58	8.00	8.68
2	8.47	8.52	8.34	8.49	8.53	8.23	8.64	7.97	8.94
3	8.54	8.38	8.31	8.50	8.32	8.20	8.37	8.02	8.66
4	8.44	8.39	8.29	8.67	8.60	8.15	8.82	7.99	8.80
5	8.38	8.41	8.26	8.57	8.30	8.19	8.50	7.93	9.03
6	8.53	8.40	8.26	8.60	8.39	8.18	8.72	8.00	8.68
7	8.54	8.40	8.29	8.48	8.50	8.11	8.58	7.98	9.00
DUP 1			8.48			8.18			8.80

MEAN	8.49	8.41	8.32	8.53	8.43	8.17	8.60	7.98	8.82
STD	0.06	0.05	0.07	0.08	0.11	0.04	0.14	0.03	0.14
minus blank	8.31	8.23	8.14	8.35	8.25	7.99	8.42	7.80	8.68

SPIKES

1/2 SAMPLE + 1/2 5 PPM KHP

	ACTUAL	ORIGINAL	EXPECTED
L1	2.73	0.37	2.68
L4	2.77	0.62	2.81
L6	2.79	0.53	2.76
L8	2.78	1.06	3.03
L9	2.73	0.50	2.75
H1	6.67	8.49	6.74
H2	6.63	8.41	6.70
H5	6.67	8.43	6.72
H7	6.61	8.60	6.80

SCALING FACTOR	0.04637 ugC/mV
ORGANIC BLANK	8.96 mV
INORGANIC BLANK	1.27 mV

VIMS
Shimadzu TOC 500, ASI 502

VIMS DOC COMPARISON STUDY
in mg/L

	L1	L2	L3	L4	L5	L6	L7	L8	L9
REP 1	0.74	0.90	1.26	1.21	0.71	0.82	0.77	1.63	0.66
2	0.70	0.83	1.23	0.95	0.72	1.00	0.44	1.45	0.78
3	0.74	0.74	1.35	0.87	0.74	0.85	0.44	1.63	0.64
4	0.68	0.74	1.20	0.87	0.78	0.87	0.48	1.63	0.68
5	0.68	0.70	1.36	0.83	0.79	0.87	0.48	1.70	0.70
6	0.59	0.64	1.29	0.87	0.82	0.87	0.42	1.67	0.69
7	0.68	0.64	1.21	0.94	0.83	0.81	0.40	1.98	0.76
MEAN	0.69	0.74	1.27	0.93	0.77	0.87	0.49	1.66	0.70
STD	0.05	0.09	0.06	0.12	0.04	0.06	0.12	0.025	0.047
minus blank	0.26	0.31	0.84	0.50	0.34	0.44	0.06	1.23	0.27
DI H2O									
BLANK									
REP 1	0.39								
2	0.47								
3	0.47								
4	0.44								
5	0.45								
6	0.39								
7	0.42								
MEAN	0.43								
STD	.03								

VIMS(2)

	H1	H2	H3	H4	H5	H6	H7	H8	H9
REP 1	8.85	8.68	8.64	9.03	8.69	8.61	9.37	8.54	9.59
2	8.81	8.68	8.67	8.92	9.17	8.63	9.06	8.50	9.56
3	8.80	8.59	8.84	8.72	9.00	8.58	8.06	8.45	9.59
4	8.90	8.54	8.66	8.74	8.96	8.63	8.39	8.32	9.58
5	9.05	8.48	8.64	8.70	8.72	8.53	8.29	8.45	9.63
6	8.70	8.59	8.72	8.78	8.65	8.49	8.30	8.54	9.69
7	8.88	8.74	8.70	8.71	8.73	8.53	8.38	8.43	9.67
MEAN	8.86	8.59	8.70	8.80	8.84	8.57	9.26	8.46	9.62
STD	0.10	0.07	0.06	0.12	0.18	0.05	0.130	0.07	0.05
minus blank	8.43	8.16	8.27	8.37	8.41	8.14	8.83	8.03	9.19

SPIKES

	sample	sample + spike	spike (3 mg/L)
L1	0.69	3.65	2.96
L2	0.74	3.98	3.24
L3	1.27	4.29	3.02
L4	0.93	4.02	3.09
L5	0.77	3.80	3.03
L6	0.87	3.95	3.08
L7	0.49	3.55	3.06
L8	0.70	3.78	3.08

Standard Reference Material True Value

DI Blank	0.04	0.38	0.14
DI SRM	3.11	3.49	3.29
SRM-Blank	3.07	3.11	3.15
Salt Blank	0.13	0.69	0.76
Salt SRM	3.19	3.98	3.83
SRM-Blank	3.06	3.29	3.07

Appendix D

DOC Samples List

Instruction for Analyses

DOC SAMPLES LIST

1. NINE SAMPLE CONTAINERS LABELED:

H1 H2 H3 H4 H5 H6 H7 H8 H9

2. NINE SAMPLE CONTAINERS LABELED:

L1 L2 L3 L4 L5 L6 L7 L8 L9

3. DI WATER BLANK WITH SAME PRESERVATION AS ABOVE SAMPLES, LABELED:
BLANK DI H2O

INSTRUCTIONS

Analyzing samples

Run seven replicates on each sample. Use the usual quality control for accuracy. If you have any standard reference material from EPA, include this in the analysis. There is no need to run additional duplicates. Please give information of which chemicals are used for the spike and SRM.

Documentation

Send a copy of the instrument standard operating procedures. Detailed information of calibration procedure for instrument is needed. Include the instrument methodology, model number, when purchased, gases used by instrument, or for sparging, and any reagents used by instrument, or in samples before loading instrument.

Data

1. Report the data for all seven replicates
2. Report any quality control results.
3. Report the standards used; what the matrix used for dilution (ie, DI water).
4. Include the regression, slope and intercept for the standards and what the chemical that was used for the standards.
5. Please send a hard copy of your results and other information. Also if it is possible to send a floppy disk with the data in an ASCII file; it will facilitate data Handling.

Appendix E

Results and graphs from Coordinated Split Sampling Program for 1990-1993 Chesapeake Bay Program

Table with Median TOC results for CSSP 1990-1991

Table with Median DOC results from CSSP 1990-1991

Figure 11.

Figure 41.

Figure 43.

Median DOC values of split sample results for 1992-1993

1992 Four-Way splits for DOC graph

1993 Four-Way splits for DOC graph

TF5.5 split sample medians with Friedman analysis results, 1990-1991
data. Data from 1991 were analyzed separately where applicable.

Parameter ¹	N ²	Laboratory Medians (mg/l)				Friedman results ³	
		DCLS	HRSD	-ODU	VIMS	X ²	P
NH4	8	0.0750	0.0950 A	0.0728 B	0.0760	18.2	<0.001
NH4 ⁴	5	0.0800	0.1000	0.0739	0.0800	6.1	<0.20
NO2	3	0.050	0.080 A	0.046 B	0.052	15.6	<0.01
NO2	2	0.035	0.055	0.027	0.032	12.7	<0.01
NO23	8	0.365	0.385	0.374	0.385	1.1	>0.70
NO23	5	0.400	0.420	0.421	0.410	7.8	<0.10
TDN	4	-	0.623 B	0.737	0.878 A	11.2	<0.01
PN	3	-	0.450	0.258	0.315	1.4	>0.30
TN	5	0.800	0.710 B	0.799 B	0.933 A	26.8	<0.001
TN	4	1.150	1.070 B	1.072 B	1.158 A	21.1	<0.001
PO4F	7	0.020 C	0.030 AB	0.037 A	0.025 BC	39.4	<0.001
PO4F	4	0.020 B	0.028	0.036 A	0.025 B	30.0	<0.001
TDP	7	0.030 B	0.060 A	0.036	0.033	24.5	<0.001
TDP	4	0.030	0.056	0.032	0.033	5.3	<0.20
PHOSP	7	0.090 A	0.090	0.078	0.064 B	15.1	<0.01
PHOSP	4	0.110 A	0.101	0.097	0.083 B	15.9	<0.01
TP	7	0.160	0.169	0.127	0.136	8.7	<0.05
TP	4	0.165	0.169	0.138	0.138	13.1	<0.01
TOC	5	3.87 B	-	7.18 A	7.43 A	21.3	<0.001

Parameter	N ¹	Laboratory Means (mg/l except CHLA & PHEA)					P value ²
		CBL	ODU	VIMS	MDHMH	DCLS	
PHOSP	7	0.0162	0.0176	0.0179	.	.	<0.30
PHOSP	7	0.0162	0.0176	0.0179	0.0175	.	<0.05
PHOSP	4	0.0157	0.0158	0.0173	0.0117	0.0192	<0.02
TP	7	0.0261	0.0267	0.0292	.	.	<0.30
TP	7	0.0261	0.0267	0.0292	0.0329	.	<0.50
TP	4	0.0272	0.0245	0.0296	0.0348	0.0375	<u>0.01</u>
DOC	6	2.8194 B	3.0483 B	3.5267 A	.	.	<u><0.001</u>
PC	8	1.3611 A	1.0688 B	1.0743 B	.	.	<u><0.001</u>
TOC	6	4.2482	4.1644 B	4.6867 A	.	.	<u><0.001</u>
TOC	3	4.1767	4.5740	3.5067	.	3.7822	<0.02
TSS	7	5.0738 B	12.0619 A	12.7905 A	.	.	<u><0.001</u>
TSS	7	5.0738 B	12.0619 A	12.7905 A	10.1429 A	.	<u><0.001</u>
TSS	4	5.2875 B	7.6833	7.7750	10.5000 A	12.0000	<u><0.001</u>
CHLA ⁴	6	.	10.6550	12.2033	12.3067	.	<0.10
PHEA ⁴	4	.	0.7767	1.7842	1.1133	.	<0.10
SI	7	0.5610 B	0.6705 A	0.6665 A	.	.	<u><0.001</u>
SI	7	0.5610 B	0.6705 A	0.6665 A	0.7048 A	.	<u><0.001</u>
SI	3	0.7789 B	0.9580 A	0.8830	0.9889 A	0.8359	<u><0.001</u>

¹ Number of cruises (sample dates) with complete data.

² Underlined values were statistically significant (P < 0.01), based on Friedman two-way ANOVA using three subsamples per cruise. Laboratory means with different letters below them also had statistically significant pairwise differences (A > B, P < 0.01).

³ Too many values were below the method detection limit to make a comparison.

⁴ Units are ug/l, not mg/l, for CHLA and PHEA.

FIGURE 11. Split sample data for Dissolved Organic Carbon (DOC), from samples collected at CB5.3 or CB4.4 (Mainstem), showing cruise means with precision bars.

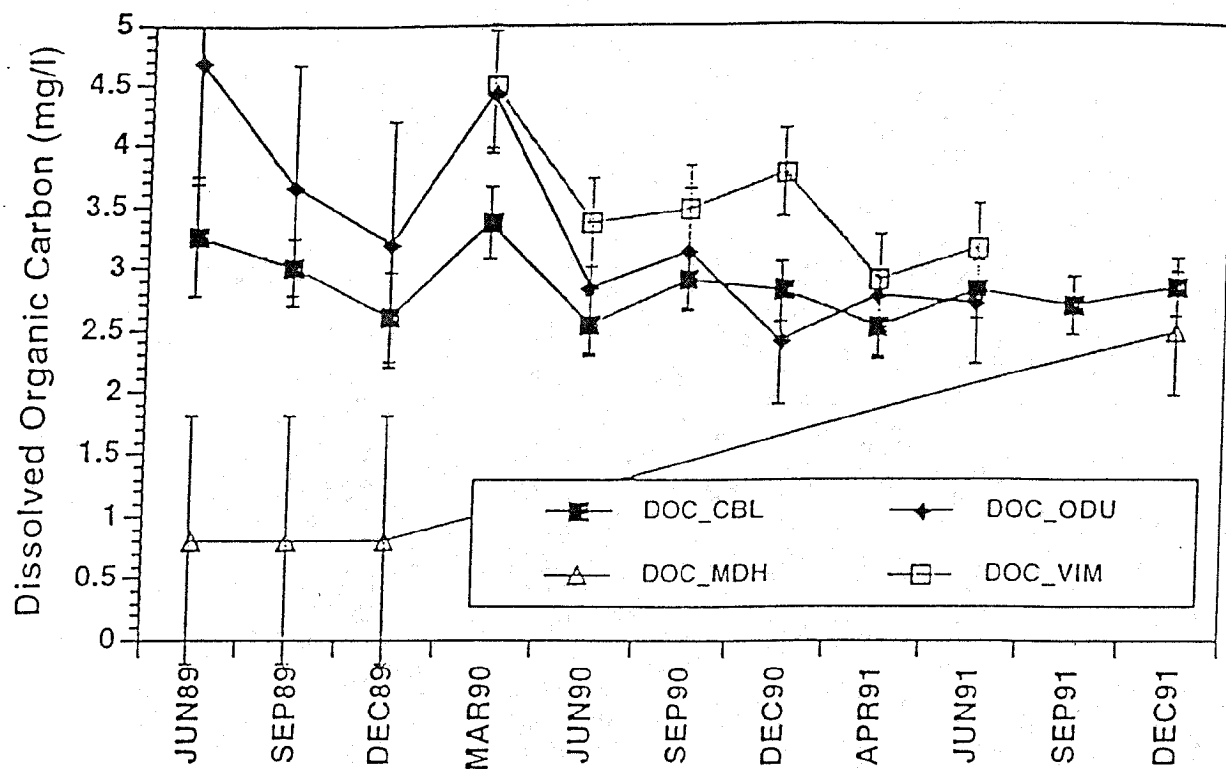


FIGURE 41. Split sample data for dissolved organic carbon (DOC), from Virginia samples collected at TF5.5 showing medians for each sample date with precision bars.



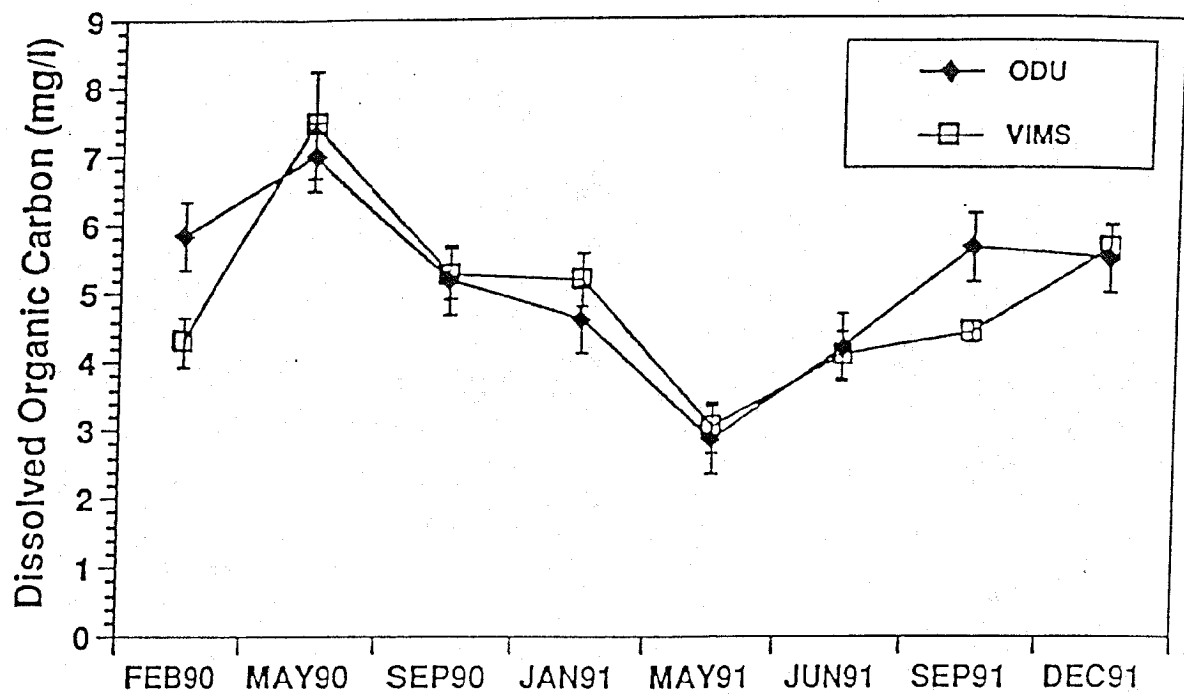


FIGURE 43. Split sample data for total organic carbon (TOC), from Virginia samples collected at TF5.5, showing medians for each sample date with precision bars.

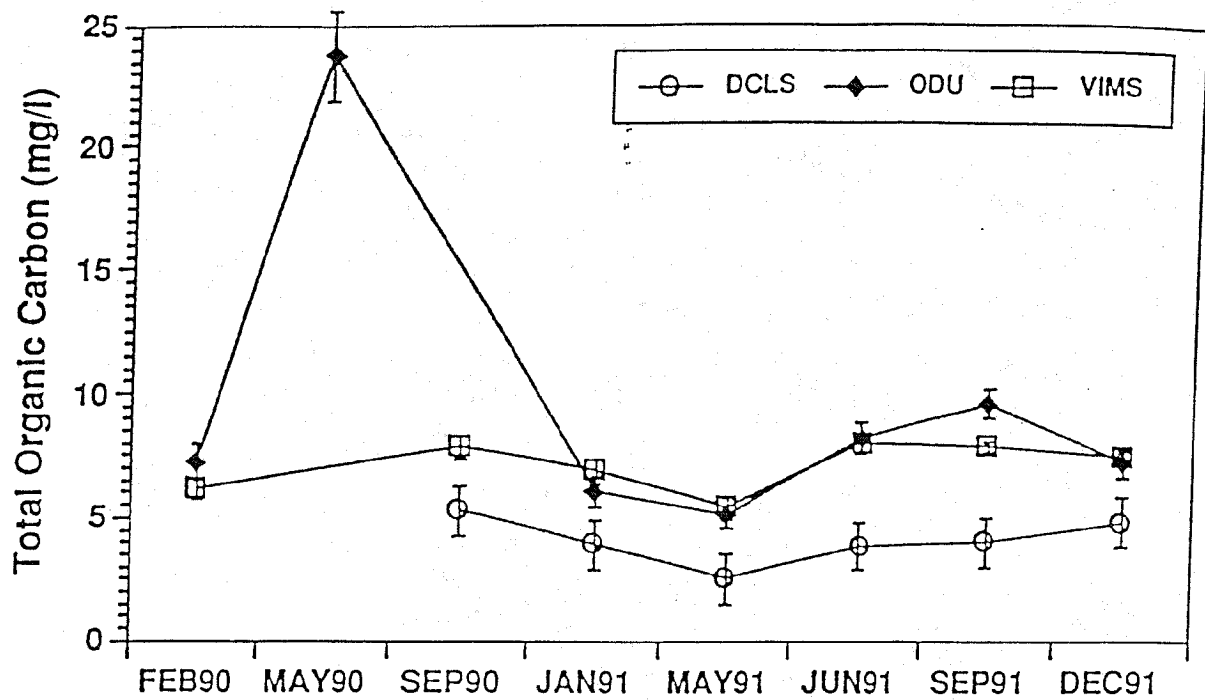

	ODU	CBL	VIMS	MDHMH	
Apr-92	3.31	2.68	3.11	1.74	
Jun-92	2.2	2.62	3.38	0.97	
Sep-92	3.015	2.54	3.79	1.03	
Dec-92	2.27	2.27	3.36	1.02	
	ODU	CBL	VIMS	MDHMH	
Feb-93	2.63	2.26	3.22	1.62	
May-93	2.3	2.74	4.1	1.85	
Aug-93	2.66	2.72	3.7	2.03	
Nov-93	3.32	2.4	3.7	2	

1992 FOUR-WAY SPLITS FOR DOC

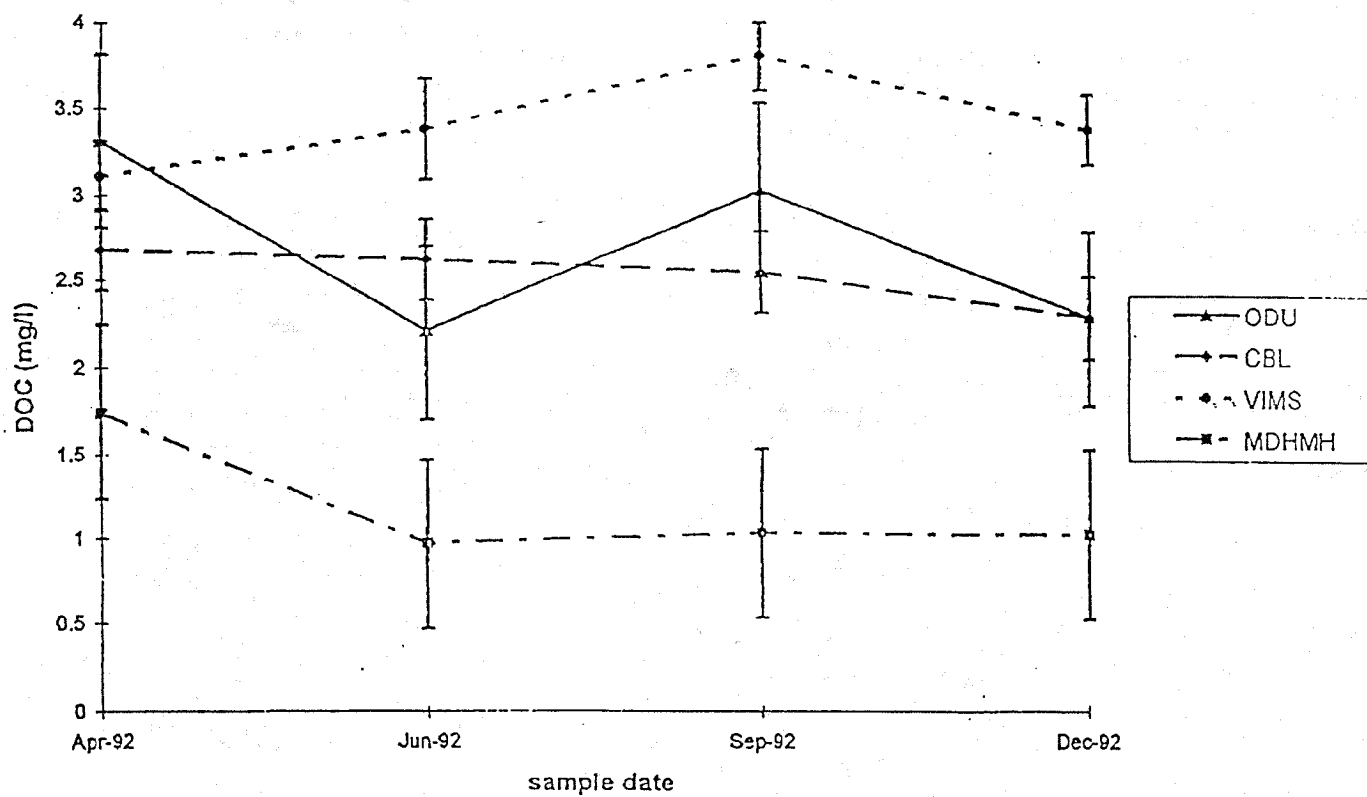

1993 FOUR-WAY SPLITS FOR DOC

