

2-45 REV 12 79



E. I. DU PONT DE NEMOURS & COMPANY
INCORPORATED
WILMINGTON, DELAWARE 19898

POLYMER PRODUCTS DEPARTMENT
EXPERIMENTAL STATION

PERSONAL AND CONFIDENTIAL

cc: A. J. Dahl - 353
B. W. Karrh - N11400
L. J. Papa - 269
Pral File
I.C.

AR226-1573

May 15, 1981

DR. B. RAMIREZ
CR&D
HASKELL LABORATORY

ANALYSIS OF BLOOD SAMPLES FOR PERFLUOROOCTANOATE
(Job No. 810-578; PRAL Nos. 81-01544 to 01552; Notebook Nos. E22514, E26238)

As requested in your letter of 4/20/81 to L. J. Papa, the nine blood samples submitted then have been analyzed for perfluorooctanoate (C₈). Results and sample identification are given in the attached table.

As noted there, the analyses were done using a gas chromatographic method specific for C₈ (Lab Method Number ES-567) but results have been reported as ppm F for comparison with total organic fluorine analyses. Precision is $\pm 10\%$ relative standard deviation over most of the concentration range, somewhat less at the lowest values. The lower limit for quantitation is 0.007 ppm F (0.01 ppm perfluorooctanoic acid), with a detection limit of ~ 0.004 ppm which can be distinguished from the reagent background but not well quantitated.

Please contact me (772-4440) or L. J. Papa (772-2745) if you have any questions regarding the analyses. General questions on blood sampling can be directed to J. W. Raines or L. F. Percival.


S. S. Stafford

Attachment
jah

Key Words:
Perfluorooctanoic Acid
Perfluorooctanoate
Blood Analysis
GC

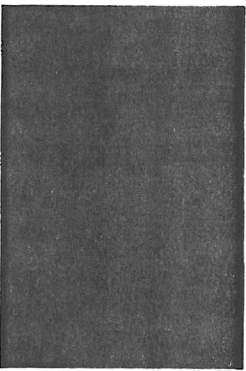
There's a world of things we're doing something about

EXP000007
EID713821

000070

TABLE I

CONCENTRATION OF PERFLUOROOCTANOATE IN BLOOD (a)

<u>Sample</u>				<u>GC Analysis</u> (b)	
<u>PRAL No.</u>	<u>Date Sampled</u>	<u>P.R.No.</u>	<u>Name</u>	<u>Date Analyzed</u>	<u>[C₈], μg F/g blood</u>
81-1544	4/20/81	-		4/21/81	0.030
81-1545	4/20/81	-		4/21/81	n.d.
81-1546	4/20/81	-		4/21/81	n.d.
81-1547	4/20/81	-		4/21/81	0.012
81-1548	4/20/81	-		4/21/81	n.d.
81-1549	4/20/81	-		4/21/81	n.d.
81-1550	4/20/81	-		4/21/81	0.0098
81-1551	4/20/81	-		4/21/81	n.d.
81-1552	4/20/81	-		4/21/81	0.0096

- (a) Analysis as described in Lab Method ES-567 ("Determination of Perfluorooctanoic Acid in Blood, Gas Chromatographic Method", S. Stafford, 4/3/81), using the packed column GC analysis with perfluoro-n-octanoic acid as calibration standard.
- (b) Although the analysis is specifically for perfluorooctanoate (acid or salts), concentrations are given in ppm fluorine for comparison with the results of total organic fluorine analyses. ($\text{ppm F} = 0.688 \times \text{ppm perfluorooctanoic acid}$) Estimated uncertainty is $\pm 10\%$ relative standard deviation. The lower limit for quantitation is $0.007 \mu\text{gF/g}$. The detection limit is $\sim 0.004 \mu\text{gF/g}$, but concentrations in that range cannot be well quantitated and are reported as < 0.007 . None detected (n.d.) is reported for samples with $[\text{C}_8] \lesssim 0.004 \text{ ppm}$, which cannot be distinguished from reagent background.

REDACTED

EXP000008
EID713822

000071