



E. I. DU PONT DE NEMOURS & COMPANY
INCORPORATED

WILMINGTON, DELAWARE 19898

POLYMER PRODUCTS DEPARTMENT
EXPERIMENTAL STATION

cc: A. J. Dahl - 353
B. W. Karrh - N11400
L. J. Papa - 269
Pral File
I.C.

Complainant's
Exhibit No. 46

PERSONAL AND CONFIDENTIAL

May 20, 1981

TO: DR. J. M. HEDGES - PPD, Circleville, Ohio

FROM: S. S. STAFFORD *S. S. Stafford*

ANALYSIS OF BLOOD SAMPLES FOR PERFLUOROOCTANOATE
(Job No. 810-633; PRAL Nos. 81-1936-1950; Notebook Nos. E22514, E26238)

As requested in your letter of 5/1/81 to L. J. Papa, the 16 blood samples submitted then have been analyzed for perfluorooctanoate (C₈). Results and sample identification are given in the attached table.

As noted there, the analyses were done using a gas chromatographic method specific for C₈ (Lab Method Number ES-567) but results have been reported as ppm F for comparison with total organic fluorine analyses. Precision is $\pm 10\%$ relative standard deviation over most of the concentration range, somewhat less at the lowest values. The lower limit for quantitation is 0.007 ppm F (0.01 ppm perfluorooctanoic acid), with a detection limit of ~ 0.004 ppm which can be distinguished from the reagent background but not well quantitated.

Please contact me (772-4440) or L. J. Papa (772-2745) if you have any questions regarding the analyses. General questions on blood sampling can be directed to J. W. Raines or L. F. Percival.

Attachment
jah

Key Words:

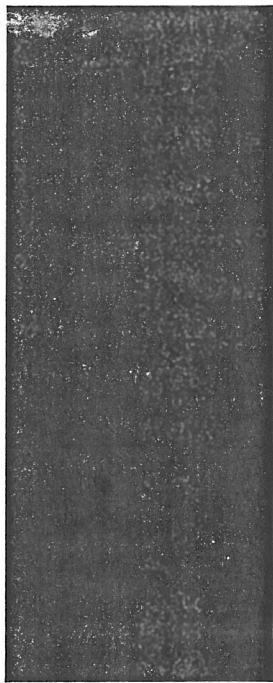
Perfluorooctanoic Acid
Perfluorooctanoate
Blood Analysis
GC

There's a world of things we're doing something about

EXP000022
EID713836

TABLE I

CONCENTRATION OF PERFLUOROOCTANOATE IN BLOOD (a)

| Sample | | | | GC Analysis | |
|----------|--------------|---------|--|---------------|---|
| PRAL No. | Date Sampled | P.R.No. | Name | Date Analyzed | [C ₈], $\mu\text{g F/g bl}$ |
| 81-1936 | 5/1/81 | 1485 |  | 5/12/81 | 0.027 |
| 81-1937 | 5/1/81 | 3518 | | 5/12/81 | n.d. |
| 81-1938 | 5/1/81 | 2430 | | 5/12/81 | 0.017 |
| 81-1939 | 5/1/81 | 3032 | | 5/12/81 | 0.015 |
| 81-1940 | 5/1/81 | 5656 | | 5/12/81 | n.d. |
| 81-1941 | 5/1/81 | 6828 | | 5/12/81 | n.d. |
| 81-1942 | 5/5/81 | 3322 | | 5/12/81 | 0.030 |
| 81-1943 | 5/4/81 | 0275 | | 5/12/81 | 0.027 |
| 81-1944 | 5/4/81 | 2718 | | 5/12/81 | 0.025 |
| 81-1945 | 5/4/81 | 3766 | | 5/12/81 | n.d. |
| 81-1946 | 5/4/81 | 4748 | | 5/13/81 | 0.017 |
| 81-1947 | 5/4/81 | 5646 | | 5/13/81 | 0.026 |
| 81-1948 | 5/4/81 | 7350 | | 5/13/81 | 0.028 |
| 81-1949 | 5/5/81 | 1505 | | 5/13/81 | 0.008 |
| 81-1950 | 5/5/81 | 5645 | | 5/13/81 | 0.053 |
| 81-1951 | 5/5/81 | 7452 | | 5/13/81 | 0.017 |

- (a) Analysis as described in Lab Method ES-567 ("Determination of perfluorooctanoic Acid in Blood, Gas Chromatographic Method", S. Stafford, 4/3/81), using the packed column GC analysis with perfluoro-n-octanoic acid as calibration standard.
- (b) Although the analysis is specifically for perfluorooctanoate (acid or salts), concentrations are given in ppm fluorine for comparison with the results of total organic fluorine analyses. ($\text{ppm F} = 0.688 \times \text{ppm perfluorooctanoic acid}$) Estimated uncertainty is $\pm 10\%$ relative standard deviation. The lower limit for quantitation is $0.007 \mu\text{gF/g}$. The detection limit is $\sim 0.004 \mu\text{gF/g}$, but concentrations in that range cannot be well quantitated and are reported as < 0.007 . None detected (n.d.) is reported for samples with $[\text{C}_8] \lesssim 0.004 \text{ ppm}$, which cannot be distinguished from reagent background.

EXP000023

REDACTED

EID713837