



3058 Research Drive  
State College, Pennsylvania 16801 USA  
Telephone: 814.272.1039  
Fax: 814.272.1019

---

## **Analytical Report**

**PFOA and PFOS Analysis of Deer Serum Samples by LC/MS/MS**

**MPI Report No. L0019345**

**Revised Report Date: 12/17/09**

---

### ***Testing Laboratory***

MPI Research, Inc.  
3058 Research Drive  
State College, PA 16801

---

### ***Requester/Project Manager***

Dena Haverland  
Dalton Utilities  
PO BOX 869  
Dalton, GA 30722  
Phone: 706-529-1010

2010 JAN -6 P 1:44

## 1 Introduction

Results are reported for the analysis of PFOS in the 3.5 yr male deer serum sample received at MPI Research from Dalton Utilities. The MPI Research study number assigned to the project is L0019345. Table I lists the target analytes quantitated for the samples.

Table I. Target Analytes for Quantitation

Compound Name	Acronym
Perfluorooctanesulfonate	C8 Sulfonate or PFOS

Note: PFOA results for both the 0.5 yr female and 3.5 yr male are reported in the original report signed on 11/19/09, as well as the PFOS results for the 0.5yr female.

## 2 Sample Receipt

Two samples were received from Dena Haverland at Dalton Utilities for this study. The samples were collected on October 02, 2009. The samples arrived on October 06, 2009 via Fedex and were logged in under MPI Research login number L0019345. The shipment was received frozen on dry ice. The samples were stored frozen at approximately -80°C from receipt until analysis. Chain-of-custody information is presented in Attachment A.

## 3 Methods - Analytical and Preparatory

### 3.1 Serum Sample Preparation

- 3.1.1. Measure 1 mL of serum sample into a 50 mL disposable centrifuge tube and fortify, if appropriate. Add 20 µL of a 50000 ng/mL WIS for a final concentration of 0.5 ng/mL.

Note: The internal standard was spiked at a higher level to allow for post extraction dilutions to be performed.

- 3.1.2 Add water to sample for a final volume of 20 mL. Cap tightly and vortex for ~1 minute.
- 3.1.3 Transfer 1 mL of the sample using a disposable pipette into 15 mL disposable centrifuge tubes. Add 5 mL of ACN and shake for ~20 minutes on a wrist action shaker.
- 3.1.4 Centrifuge tubes at ~3000 rpm for ~5 minutes. Carefully decant supernatant into a 50 mL disposable centrifuge tube and add 35 mL of water.
- 3.1.5 Place the unconditioned SPE columns on the vacuum manifold. Condition the SPE columns by passing ~10 mL of methanol through the column followed by ~5 mL of water. The washes may be pulled through the SPE column using vacuum at a flow rate of ~1 drop/sec or may be allowed to pass through the column unaided. Discard all washes. Do not allow the column to dry.
- 3.1.6 Load the sample onto a conditioned SPE column. Discard the eluate. Any analyte residues will be trapped on the SPE column at this point.

3.1.7 Elute with 2 mL of methanol. Collect 2 mL of elute into a graduated 15 mL centrifuge tube.

Note: Post extraction dilutions were prepared in methanol.

### 3.2 Sample Analysis by LC/MS/MS

In High Pressure Liquid Chromatography (HPLC), an aliquot of extract is injected and passed through a liquid-phase chromatographic column. Based on the affinity of the analyte for the stationary phase in the column relative to the liquid mobile phase, the analyte is retained for a characteristic amount of time. Following HPLC separation, mass spectrometry provides a rapid and accurate means for analyzing a wide range of organic compounds. Molecules are ionized, fragmented, and detected. The ions characteristic of the compounds are observed and quantitated against external calibration standards.

An HP1100 system interfaced to an Applied Biosystems API 4000 LC/MS/MS was used to analyze the sample extracts for quantitation. A gradient elution through a Phenomenex Luna 3 $\mu$  C8(2) Mercury, 20 x 4.0 mm column was used for separation.

The following gradient was performed:

Mobile Phase (A): 2mM Ammonium Acetate in Water  
Mobile Phase (B): Methanol

Time	%A	%B
0.0	90	10
0.5	90	10
2.0	10	90
5.0	10	90
5.1	0	100
6.0	0	100
6.1	90	10
10.0	90	10

The following parameters were used for operation of the mass spectrometer:

Parameter	Setting
Ionization Mode	Electrospray
Polarity	Negative
Transitions Monitored	499→80 (PFOS) 503→80 (Internal Std. <sup>13</sup> C PFOS (m+4))
Gas Temperature	450°C

## 4 Analysis by LCMSMS

### 4.1 Calibration

For the serum sample analysis, a 6-point calibration curve was analyzed throughout the analytical sequence for PFOS. The calibration points were prepared at 0.1, 0.2, 0.5, 1.0, 2.0, 5.0 ng/mL (ppb) containing 1.0 ng/mL <sup>13</sup>C-PFOS (m+4).

The ratio of the analyte concentration to the IS concentration versus the ratio of the analyte instrument response (area) to the IS response (area) was plotted for each point. Using linear regression with 1/x weighting, the slope, y-intercept and coefficient of determination ( $r^2$ ) were determined. A calibration curve is acceptable if  $r^2 \geq 0.985$ .

For the results reported here, calibration criteria were met. The calibration curves are included in the raw data in Attachment C.

#### **4.2 Laboratory Control Spikes**

Laboratory control spikes in the analytical set were prepared during each extraction set by adding a known concentration of the analyte to deer serum controls. Laboratory control spikes are used to assess method accuracy. The laboratory control spikes must show recoveries between 70-130% or the data is rejected. For the results reported here, the laboratory control spikes were within the acceptable range. Laboratory control spike recoveries are given in Attachment B.

#### **4.3 Matrix Spikes**

One matrix spike was prepared by adding a known concentration of the target analyte to a sample. Matrix spikes are used to assess method accuracy in the matrix. The matrix spikes should show recoveries between 70-130%. For the results reported here, the matrix spike was within the acceptable range with the exceptions of.

#### **4.4 Laboratory Duplicates**

One sample was prepared in duplicate and analyzed. Duplicate results are given along with the sample results in Attachment B.

### **5 Data Summary**

Please see Attachment B for a detailed listing of the analytical results. For the serum samples the results are reported in parts per billion (ng/mL) on an as-received basis.

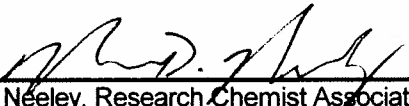
### **6 Data/Sample Retention**

Samples are disposed of 60 days after the report is issued unless otherwise specified by the project manager. All electronic data is archived on retrievable media and hard copy reports are stored in data folders maintained by MPI Research. Hardcopy data is stored for a minimum of five years. The client will be notified 30 days prior to the disposal of hardcopy data.

## 7 Attachments

- 7.1 Attachment A: Chain of Custody
- 7.2 Attachment B: Analytical Results
- 7.3 Attachment C: Raw Analytical Data for Water

## 8 Signatures

  
\_\_\_\_\_  
Mark Neeley, Research Chemist Associate II      12-17-09  
Date

  
\_\_\_\_\_  
Robert Zhu, Manager, Analytical      12/17/09  
Date



A







Mattawan (Corporate Headquarters)  
54943 North Main Street  
Mattawan, MI 49071-9399  
(269) 668-3336 Phone  
(269) 668-4151 Fax

State College  
3058 Research Drive  
State College, PA 16801  
(814) 272-1039 Phone  
(814) 231-1580 Fax

## Login

Login Group: L0019345

Login #:	19459	Conform COC Sample:	True
Project:	P0005195	Conform COC:	True
Company Name:	Dalton Utilities	Conform Sample:	True
Submitted By:	Dena Haverland	Conform Request:	True
Login Type:	Immediate Receipt of Samples		
Started:	True		
Date Start:	10/27/2009		
Due Date:	11/06/2009		
Login Initiated:	10/27/2009		
Received By:	Ammerman, Mark		
Spread Sample:			
Label:			
MPI SD/PI:	Zhu, Xiang		
Project Title/Type:	PFOA and PFOS Analysis of Serum Samples By LC/MS/MS / ROUTINE		
Login Notes:			

## Packages / Containers

Package	Carton	Date / Condition	Shipper / ID	Temp. Control/Temp.	Direction / Handled By
K0022041		Received Date: 10/6/09 10:25 Package & Contents Uncompromised	FEDEX 8694 2057 8178	Dry Ice -79.2	RECEIVED Ammerman, Mark

Container #	Gross Weight	pH	Container Type	Preservative	Mfg. Lot	Mfg. ID
604	3.10 g		2 ml clear plst vial	NONE		
C0457605	3.10 g		2 ml clear plst vial	NONE		
C0457606	3.20 g		2 ml clear plst vial	NONE		
C0457607	3.10 g		2 ml clear plst vial	NONE		
C0457608	3.10 g		2 ml clear plst vial	NONE		
C0457609	3.10 g		2 ml clear plst vial	NONE		
C0457610	3.20 g		2 ml clear plst vial	NONE		
C0457611	3.40 g		2 ml clear plst vial	NONE		
C0457612	3.10 g		2 ml clear plst vial	NONE		
C0457613	3.10 g		2 ml clear plst vial	NONE		
C0457614	3.10 g		2 ml clear plst vial	NONE		
C0457615	3.20 g		2 ml clear plst vial	NONE		
C0457616	3.10 g		2 ml clear plst vial	NONE		
C0457617	3.10 g		2 ml clear plst vial	NONE		
C0457618	3.10 g		2 ml clear plst vial	NONE		
C0457619	3.10 g		2 ml clear plst vial	NONE		
C0457620	3.10 g		2 ml clear plst vial	NONE		
C0457621	3.10 g		2 ml clear plst vial	NONE		
C0457622	3.20 g		2 ml clear plst vial	NONE		
C0457623	3.10 g		2 ml clear plst vial	NONE		



## Samples

<u>Sample ID</u>	<u>Container</u>	<u>Matrix</u>	<u>System</u>	<u>System Matrix</u>	<u>Sample</u>	<u>Date Sampled</u>	<u>Date Due</u>
L0019345-0001		LIQUID	Deer	Serum	Deer #6 0.5 yr female-serum	10/02/2009	11/06/2009
	C0457613						
	C0457612						
	C0457611						
	C0457610						
	C0457609						
	C0457608						
	C0457607						
	C0457606						
	C0457605						
	C0457604						
L0019345-0002		LIQUID	Deer	Serum	Deer #7 3.5 yr male-serum	10/02/2009	11/06/2009
	C0457617						
	C0457623						
	C0457622						
	C0457621						
	C0457620						
	C0457619						
	C0457618						
	C0457616						
	C0457615						
	C0457614						

Login Reviewed By: \_\_\_\_\_

Date/Time: \_\_\_\_\_



# MPI RESEARCH

MPI Research Contact: Daniel Wright

## Send Report To:

Company: Dalton Utilities  
 Address: 1200 VD Parrott JK Parkway, PO Box 869  
 City, State, ZIP: Dalton, GA 30722-0869  
 Attention: Dena Haverland  
 Phone #: 706-529-1010  
 Fax #: 706-529-1271  
 Email: dhaverland@dutil.com  
 Study/Job #: \_\_\_\_\_  
 Signature/Date: \_\_\_\_\_  
 Printed Name: \_\_\_\_\_

## Sample Submittal

Please fax this form before sending samples.

Please send samples to shipping and receiving:  
 3048 Research Drive, State College, PA 16801  
 T: (814) 272-1039 • F: (814) 272-1019

## Turnaround time (TAT) requirements:

Results Due Date: 30 days  
 Preliminary Results Format: Verbal ☒ Email ☐ Fax  
 Report Due Date: 30 days

## Storage Conditions

Room temperature  
 Refrigerator  
☒ Freezer  
 Ultra Low freezer  
 Desiccated  
 Lighting required

Stability (°C/%RH): \_\_\_\_\_

Stability time period: \_\_\_\_\_

## Safety Information

Special handling: \_\_\_\_\_  
 MSDS attached  
 Controlled substance: \_\_\_\_\_  
 HAZARDS: \_\_\_\_\_  
 Please fill in the diamond HMIS/NFPA  
 (0 4) if appropriate

	Client ID# Description	Lot/ Control #	Amt. Sent/ Weight	# of Bottles	Matrix	Date & Time	Tests Requested
1	Deer #6 0.5 yr female-Serum		10ml	10	deer	10/2/09 1:08AM	PFOA/PFOS
2	Deer #6 0.5 yr female-Muscle		as requested	1 bag	deer	10/2/09 2:28AM	PFOA/PFOS
3	Deer #6 0.5 yr female-Liver		Whole	1 bag	deer	10/2/09 2:30AM	PFOA/PFOS
4	Deer #7 3.5 yr Male-Serum		10ml	10	deer	10/2/09 1:45AM	PFOA/PFOS
5	Deer #7 3.5 yr Male-Muscle		as requested	1 bag	deer	10/2/09 2:45AM	PFOA/PFOS
6	Deer #7 3.5 yr Male-Liver		Whole	1 bag	deer	10/2/09 2:48AM	PFOA/PFOS
7							
8							
9							
10							

PO #: \_\_\_\_\_

Relinquished by	Date	Time	Received by	Date	Time
<u>Daniel Karam</u>	<u>10/5/09</u>	<u>6:30 PM</u>	<u>[Signature]</u>	<u>10/6/09</u>	<u>1025</u>

Notes:

\_\_\_\_\_  
 \_\_\_\_\_  
 \_\_\_\_\_  
 \_\_\_\_\_  
 \_\_\_\_\_



## TEMPORARY SAMPLE STORAGE FORM

To be completed during ExyLIMS Login

Project #: P5195

Login #: L19345

Initials / Date: MA 10/27/09

One form to be completed for each package

Date / Time Received: 10/06/09 1025

Received By: Mark Ammer

Shipper: FedEx

Shipper Package ID: 8694 2057 8178

Temperature (deg C) / Thermometer ID: -79.2 / D0001775

Temperature Control Method: dry ice-active

Temporary Storage Location: freezer 32

Condition of sample(s):

- ☒ Good – Package and contents uncompromised  
☐ Fair – Package damaged / contents uncompromised  
☐ Poor – Package and contents compromised

Notes:





B







**MPI**  
RESEARCH

3058 Research Drive  
State College, Pennsylvania 16801 USA  
Telephone: 814.272.1039  
Fax: 814.272.1019

## Analytical Report

### Summary of Fluorochemical Residues in Serum Samples

Sample ID	PFOS Perfluorooctanesulfonate
	Analyte Found (ng/mL, ppb)
Deer # 7 3.5 yr male-serum	1140
Deer # 7 3.5 yr male-serum*	1120

\*Laboratory Duplicate

ND = Not detected = Response is below the LOD of 1.0 ng/mL (ppb).

NQ = Not quantifiable = Response is between the LOD and the LOQ of 10 ng/mL (ppb).



3058 Research Drive  
State College, Pennsylvania 16801 USA  
Telephone: 814.272.1039  
Fax: 814.272.1019

## Recovery Summary of Fluorochemical Residues in Serum Samples

Sample Description	Amount Spiked (ng/mL)	Amt Found in Sample (ng/mL)	PFOS	
			Amount Recovered (ng/mL)	Recovery (%)
LCS A (Data set 120709A) 1000 ng/mL	1000	ND	1050	105
LCS B (Data set 120709A) 1000 ng/mL	1000	ND	1100	110
Deer # 7 3.5 yr male-serum (L19345-2 Spk C, 1000 ng/mL Lab Spike)	1000	1140	2130	99

ND = Not detected = Response is below the LOD of 1.0 ng/mL.

NQ = Not quantifiable = Response is between the LOD and the LOQ of 10 ng/mL.