
Analytical Report

Fluorochemical Characterization of Aqueous Samples

MPI Report No. L0019716

Testing Laboratory

MPI Research, Inc.
3058 Research Drive
State College, PA 16801

Requester

Jason Governo
University of Georgia
Driftmier Engineering Center
Athens, GA 30602
Phone: 678-794-6664

1 Introduction

Results are reported for the analysis of eleven water samples received at MPI Research from the University of Georgia. The MPI Research study number assigned to the project is L0019716. Table I lists the target analyte quantitated for the sample.

Table I. Target Analyte for Quantitation

Compound Name	Acronym
Perfluorobutyric Acid	C4 Acid
Perfluoropentanoic Acid	C5 Acid
Perfluorohexanoic Acid	C6 Acid
Perfluoroheptanoic Acid	C7 Acid
Perfluorooctanoic Acid	C8 Acid
Perfluorononanoic Acid	C9 Acid
Perfluorodecanoic Acid	C10 Acid
Perfluoroundecanoic Acid	C11 Acid
Perfluorododecanoic Acid	C12 Acid
Perfluorotridecanoic Acid	C13 Acid
Perfluorotetradecanoic Acid	C14 Acid
Perfluorobutanesulfonate	C4 Sulfonate
Perfluorohexanesulfonate	C6 Sulfonate
Perfluoroheptanesulfonate	C7 Sulfonate
Perfluorooctanesulfonate	C8 Sulfonate
Perfluorodecanesulfonate	C10 Sulfonate
Perfluorooctanesulfonamide	FOSA
2-(N-methylperfluorooctanesulfonamido)	MeFOSAA
2-(N-ethylperfluorooctanesulfonamido)	EtFOSAA
N-methylperfluorooctane	MeFOSE
N-ethylperfluorooctane	EtFOSE
6:2 Fluorotelomer alcohol	6:2 FTOH
7:2 sFluorotelomer alcohol	7:2s FTOH
8:2 Fluorotelomer alcohol	8:2 FTOH
10:2 Fluorotelomer alcohol	10:2 FTOH

2 Sample Receipt

Eleven water samples were received from the client cooled with wet ice on 12/08/09 and given the MPI Research login number of L0019716. The samples were stored refrigerated from receipt until analysis. Chain-of-custody information is presented in Attachment A.

3 Methods - Analytical and Preparatory

3.1 Water Sample Preparation for LC/MS/MS

Ten milliliters of sample was transferred into a 50 mL centrifuge tube. Ten milliliters of acetonitrile was added to the sample. After shaking, the sample was sonicated for approximately 2 hours then centrifuged at 3000 rpm for ~10 minutes. A 1 mL portion of the supernatant was transferred to an autosampler vial and fortified with an internal standard solution. The samples were then analyzed using electrospray LC/MS/MS.

3.2 Water Sample Preparation for GC/MS

One hundred milliliters of sample was transferred into a 250 mL polypropylene bottle. Forty milliliter of methyl tert-butyl ether (MTBE) was added to the bottle. The bottle was capped and then shaken for one hour on a reciprocation shaker. The content of the bottle was poured into a 250 mL separatory funnel. The bottle was rinsed with approximately 10 mL fresh MTBE and the rinsate was added to the separatory funnel. The aqueous phase in the funnel was discarded. The organic phase was collected in a 300 mL flask and dried with sodium sulfate. The dried organic phase was then quantitatively transferred into a 50 mL polypropylene centrifuge tube and concentrated to 1 mL using a nitrogen evaporator. The extracted was transferred into a 2 mL GC vial and 10 μ L of internal standard was added. The sample was then analyzed by GC/MS.

3.3 Sample Analysis by LC/MS/MS

In High Pressure Liquid Chromatography (HPLC), an aliquot of extract is injected and passed through a liquid-phase chromatographic column. Based on the affinity of the analyte for the stationary phase in the column relative to the liquid mobile phase, the analyte is retained for a characteristic amount of time. Following HPLC separation, mass spectrometry provides a rapid and accurate means for analyzing a wide range of organic compounds. Molecules are ionized, fragmented, and detected. The ions characteristic of the compounds are observed and quantitated against calibration standards.

An HP1100 system interfaced to an Applied Biosystems API 5000 LC/MS/MS was used to analyze the sample extracts for quantitation. A gradient elution through a Phenomenex Luna 3 μ C8(2) Mercury, 20 x 4.0 mm column was used for separation.

The following gradient was performed for C4-C14 acids, PFBS, PFHS, PFHpS, PFOS, PFDS, FOSA, and ¹³C PFOA (m+4):

Mobile Phase (A):	2mM Ammonium Acetate in Water	
Mobile Phase (B):	Methanol	
<u>Time</u>	<u>%A</u>	<u>%B</u>
0.0	90	10
0.5	90	10
2.0	10	90
5.0	10	90

5.1	0	100
6.0	0	100
6.1	90	10
10.0	90	10

The following gradient was performed for MeFOSAA, EtFOSAA, MeFOSE and EtFOSE:

Mobile Phase (A): 2mM Ammonium Acetate in Water
Mobile Phase (B): Methanol

Time	%A	%B
0.0	75	25
0.5	75	25
2.0	10	90
5.0	10	90
5.1	0	100
6.0	0	100
6.1	75	25
10.0	75	25

The following parameters were used for operation of the mass spectrometer:

Parameter	Setting
Ionization Mode	Electrospray
Polarity	Negative
Transitions Monitored	213→169 (C4 Acid) 263→219 (C5 Acid) 313→269 (C6 Acid) 363→319 (C7 Acid) 413→369 (C8 Acid) 463→419 (C9 Acid) 513→469 (C10 Acid) 563→519 (C11 Acid) 613→569 (C12 Acid) 663→619 (C13 Acid) 713→669 (C14 Acid) 299→80 (PFBS) 399→80 (PFHS) 449→99 (PFHpS) 499→80 (PFOS) 599→99 (PFDS) 498→78 (FOSA) 217→172 (Internal Std. ¹³ C PFBA (m+4)) 415→370 (Internal Std. ¹³ C PFOA (m+2)) 515→470 (Internal Std. ¹³ C PFDA (m+2)) 503→80 (Internal Std. ¹³ C PFOS (m+4)) 417→372 (Surrogate ¹³ C PFOA (m+4)) 570→419 (MeFOSAA) 584→419 (EtFOSAA) 616→59 (MeFOSE) 630→59 (EtFOSE)
Gas Temperature	400°C

3.4 Sample Analysis by GC/MS

The extracts were injected into a gas chromatograph (GC) equipped with a narrow bore capillary column and mass selective detector. The GC was temperature programmed to separate the analytes, and the analytes eluted from the column were introduced to the mass selective detector and identified by comparing retention times and abundances of characteristic masses to that of known standards. Sample concentration was calculated by comparing the response of the characteristic mass relative to that of the calibration curve.

The GC/MS system was operated using the following conditions:

Instrument	Hewlett-Packard model 6890 Series Gas Chromatograph/model 5973 Mass Selective Detector
Column	HP-1, 30 m x 0.25 mm ID, 1.00 μ m df
Oven Temperature	Hold at 60°C for 4 min., ramp at 20°C/min. to 140°C, ramp at 40°C to 240°C, hold for 5 minutes
Injector Temperature	200 °C
Transfer Line Temperature	280 °C
Carrier Gas	Helium
Column Flow	1.0 mL/min, Constant
Injection Mode	Pulsed Splitless, 30psi for 1.5 min.
Injection Liner	4 mm ID Single Gooseneck packed with glass wool
Injection Purge Delay	1.5 min.
Purge Flow to Split Vent	50 mL/min.
Injection Volume	2 μ L
Electron Multiplier Voltage	From ATUNE + 306V
MS Acquisition Mode	SIM
Ions Monitored	MFOET (Internal Standard): m/z 448, m/z 466 8:1 FTOH (Surrogate): m/z 363, m/z 431 6:2 FTOH: m/z 344, m/z 363 7:2s FTOH: m/z 319, m/z 355 8:2 FTOH: m/z 405, m/z 463 10:2 FTOH: m/z 505, m/z 544
Dwell Time	50ms for each ion
MS Temperature	Quad: 150 °C, Source: 230 °C

4 Analysis by LCMSMS

4.1 Calibration

A 9-point calibration curve was analyzed at the beginning of the analytical sequence for C4-C14 acids, PFBS, PFHS, PFHpS, PFOS, PFDS, FOSA, and ¹³C PFOA (m+4). A continuing calibration verification (CCV) standard (0.250 ng/mL) was used to verify the accuracy of the calibration curve for the duration of the analytical run. At the minimum every tenth sample was a CCV, not including solvent blanks. The calibration curve and the last passing CCV (70-130%) will then bracket acceptable samples. The calibration points were prepared at 0.0125, 0.025, 0.050, 0.100, 0.250, 0.500, 1.0, 2.5 and 5.0 ng/mL (ppb) for LC/MS/MS analysis. The ratio of the analyte concentration to the IS concentration versus the ratio of the analyte instrument response (area) to the IS response (area) was plotted for each point. Using linear

regression with 1/x weighting, the slope, y-intercept and coefficient of determination (r^2) were determined. A calibration curve is acceptable if $r^2 \geq 0.985$.

A 9-point calibration curve was analyzed at the beginning of the analytical sequence for MeFOSAA, EtFOSSA, MeFOSE and EtFOSE. A continuing calibration verification (CCV) standard (0.250 ng/mL) was used to verify the accuracy of the calibration curve for the duration of the analytical run. At the minimum every tenth sample was a CCV, not including solvent blanks. The calibration curve and the last passing CCV (70-130%) will then bracket acceptable samples. The calibration points were prepared at 0.0125, 0.025, 0.050, 0.100, 0.250, 0.500, 1.0, 2.5 and 5.0 ng/mL (ppb) for LC/MS/MS analysis. The instrument response versus the concentration was plotted for each point. Using quadratic regression with 1/x weighting, the X variable 1 (a), X variable 2 (b), intercept (c) and coefficient of determination (r^2) were determined. A calibration curve is acceptable if $r^2 \geq 0.985$.

For the results reported here, calibration criteria were met. The calibration curve is included in the raw data in Attachment C.

4.2 Surrogates and Internal Standards

^{13}C labeled-perfluorooctanoic acid (^{13}C PFOA (m+4)) is used as a surrogate for the water samples.

^{13}C PFOA (m+4) recoveries can be found in Attachment B.

^{13}C PFBA (m+4) is used as the internal standard for the water samples for C4 – C6 Acids.

^{13}C PFOA (m+2) is used as the internal standard for the water samples for C7 – C9 Acids.

^{13}C PFDA (m+2) is used as the internal standard for the water samples for C10 – C14 Acids.

^{13}C PFOS (m+4) is used as the internal standard for the water samples for PFBS, PFHS, PFOS and FOSA

4.3 Laboratory Control Spikes

Laboratory control spikes in the analytical set were prepared during each extraction set by adding a known concentration of the analyte to laboratory reagents. Laboratory control spikes are used to assess method accuracy. The laboratory control spikes must show recoveries between 70-130% or the data is rejected. For the results reported here, the laboratory control spikes were within the acceptable range. Laboratory control spike recoveries are given in Attachment B.

4.4 Matrix Spikes

Two matrix spikes were prepared for the water samples, one for C4-C14 acids, PFBS, PFHS, PFHpS, PFOS, PFDS, FOSA, and ^{13}C PFOA (m+4) analysis, and one for MeFOSAA, EtFOSSA, MeFOSE and EtFOSE analysis. They were prepared by adding a known concentration of the target analyte to a separate sample. Matrix spikes are used to assess method accuracy in the matrix. The matrix spike should show recoveries between 70-130%. For the results reported here, the matrix spike was within the acceptable range with the exception of:

L19716-1 (16-1) Spk C at 1.0 ng/mL for C11 acid, C13 acid, and C14 acid, which gave high recoveries.

Since the matrix spike recoveries were bias high and the samples were non-detected the data is considered acceptable.

Matrix spike recoveries are given in Attachment B.

4.5 Duplicate

Laboratory duplicates were not performed as part of this study.

5 Analysis by GC/MS

5.1 System Suitability and Calibration

Three system suitability standards were analyzed at the beginning of the analytical sequence. The %RSD of the peak area of each analyte should be ≤ 20 .

A 6-point calibration was analyzed. The calibration standard analyses were interspersed throughout the analytical sequence. The calibration points were prepared at 0.1, 0.2, 0.5, 1.0, 2.0 and 5.0 $\mu\text{g/mL}$, which are equivalent to 1, 2, 5, 10, 20 and 50 $\mu\text{g/L}$ (ppb) in samples. A calibration curve is acceptable if $r^2 \geq 0.985$.

For the results reported here, system suitability and calibration criteria were met. The system suitability and calibration curve are included in the raw data in Attachment C.

5.2 Surrogate and Internal Standard

1H,1H-Perfluoro-1-nonanol (8:1 FTOH) was used as surrogate standard. The recoveries of 8:1 FTOH can be found in Attachment B.

2-Perfluorooctyl-[1,1-2H2]-[1,2-13C2]-ethanol (MFOET) was used as internal standard.

5.3 Laboratory Control Spikes

Laboratory control spikes in the analytical set were prepared during each extraction set by adding a known concentration of the analyte to laboratory reagents. Laboratory control spikes are used to assess method accuracy. The laboratory control spikes must show recoveries between 50-120%. For the results reported here, the laboratory control spikes were within the acceptable range. Laboratory control spike recoveries are given in Attachment B.

5.4 Duplicate

A duplicate sample was not analyzed for this sample login.

5.5 Matrix Spike

A matrix spike was not analyzed for this sample login.

6 Data Summary

6.1 GC/MS Sample Results

Five samples (16-3, 18-3, 19-3, 21-3, and 3D3) were re-extracted due to surrogate recoveries lower than 50%. Upon re-extraction the surrogate recoveries of the five samples were still below 50%. The low recovery may be due the matrix effect of the sample. For each of the five samples that were re-extracted, the extraction with the higher recovery was reported.

The results are reported in parts per billion (ng/mL) on an as-received basis.

Please see Attachment B for a detailed listing of the analytical results.


7 Data/Sample Retention

Samples are disposed of one month after the report is issued unless otherwise specified. All electronic data is archived on retrievable media and hard copy reports are stored in data folders maintained by MPI Research. Hardcopy data is stored for a minimum of five years. The client will be notified 30 days prior to the disposal of hardcopy data.


8 Attachments

- 7.1 Attachment A: Chain of Custody
- 7.2 Attachment B: Analytical Results
- 7.3 Attachment C: Raw Analytical Data (LC/MS/MS)
- 7.4 Attachment D: Raw Analytical Data (GC/MS)

9 Signatures



Mark Neeley, Research Chemist Associate II
1-22-10
Date

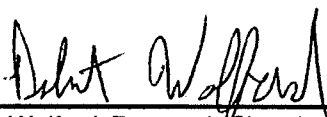

Robert Wolford, Research Chemist Associate I
1-22-10
Date



Xiang Zhu, Manager, Analytical
1/22/2010

B



Fluorochemical Residues in Water Samples By LC/MS/MS

Sample ID: L0019716-1; 16-1
Date of Extraction: 12/16/2009

Analyte	Amount Found (ng/mL)	LOD (ng/mL)	Analysis Date
C4 Acid- Perfluorobutyric Acid	NQ	0.0125	12/16/2009
C5 Acid- Perfluoropentanoic Acid	0.102	0.0125	12/16/2009
C6 Acid- Perfluorohexanoic Acid	0.238	0.0125	12/16/2009
C7 Acid- Perfluoroheptanoic Acid	0.0289	0.0125	12/16/2009
C8 Acid- Perfluorooctanoic Acid	0.0320	0.0125	12/16/2009
C9 Acid- Perfluorononanoic Acid	ND	0.0125	12/16/2009
C10 Acid- Perfluorodecanoic Acid	0.0255	0.0125	12/16/2009
C11 Acid- Perfluoroundecanoic Acid	ND	0.0125	12/16/2009
C12 Acid- Perfluorododecanoic Acid	ND	0.0125	12/16/2009
C13 Acid- Perfluorotridecanoic Acid	ND	0.0125	12/16/2009
C14 Acid- Perfluorotetradecanoic Acid	ND	0.0125	12/16/2009
PFBS- Perfluorobutanesulfonate	ND	0.0125	12/17/2009
PFHS- Perfluorohexanesulfonate	ND	0.0125	12/17/2009
PFOS- Perfluorooctanesulfonate	0.162	0.0125	12/17/2009
FOSA- Perfluorooctane sulfonamide	ND	0.0125	12/17/2009
PFHpS- Perfluoroheptanesulfonate	ND	0.0125	12/17/2009
PFDS- Perfluorodecanesulfonate	ND	0.0125	12/16/2009
MeFOSAA- 2(N-Methylperfluorooctanesulfonamido) acetic acid	ND	0.0125	12/16/2009
EtFOSAA- 2(N-Ethylperfluorooctanesulfonamido) acetic acid	ND	0.0125	12/16/2009
MeFOSE- N-Methylperfluorooctane sulfonamidoethanol	ND	0.0125	12/16/2009
EtFOSE- N-Ethylperfluorooctane sulfonamidoethanol	ND	0.0125	12/16/2009

ND = Not detected = Response is below the LOD of 0.0125 ng/mL (ppb).

NQ = Not quantifiable = Response is between the LOD and the LOQ of 0.0250 ng/mL (ppb).



3058 Research Drive
State College, Pennsylvania 16801 USA
Telephone: 814.272.1039
Fax: 814.272.1019

Analytical Report

Fluorotelomer Analysis by GC/MS

Sample ID: L0019716-3; 16-3

Date of Extraction: 1/12/2010

Date Analyzed: 1/12/2010

Analyte	Amount Found (ng/mL)	LOQ (ng/mL)
7-2s FTOH	ND	1.00
6-2 FTOH	3.40	1.00
8-2 FTOH	ND	1.00
10-2 FTOH	ND	1.00

Fluorochemical Residues in Water Samples By LC/MS/MS

Sample ID: L0019716-4; 17-1
Date of Extraction: 12/16/2009

Analyte	Amount Found (ng/mL)	LOD (ng/mL)	Analysis Date
C4 Acid- Perfluorobutyric Acid	ND	0.0125	12/16/2009
C5 Acid- Perfluoropentanoic Acid	ND	0.0125	12/16/2009
C6 Acid- Perfluorohexanoic Acid	ND	0.0125	12/16/2009
C7 Acid- Perfluoroheptanoic Acid	ND	0.0125	12/16/2009
C8 Acid- Perfluorooctanoic Acid	ND	0.0125	12/16/2009
C9 Acid- Perfluorononanoic Acid	ND	0.0125	12/16/2009
C10 Acid- Perfluorodecanoic Acid	ND	0.0125	12/16/2009
C11 Acid- Perfluoroundecanoic Acid	ND	0.0125	12/16/2009
C12 Acid- Perfluorododecanoic Acid	ND	0.0125	12/16/2009
C13 Acid- Perfluorotridecanoic Acid	ND	0.0125	12/16/2009
C14 Acid- Perfluorotetradecanoic Acid	ND	0.0125	12/16/2009
PFBS- Perfluorobutanesulfonate	ND	0.0125	12/17/2009
PFHS- Perfluorohexanesulfonate	ND	0.0125	12/17/2009
PFOS- Perfluorooctanesulfonate	ND	0.0125	12/17/2009
FOSA- Perfluorooctane sulfonamide	ND	0.0125	12/17/2009
PFHpS- Perfluoroheptanesulfonate	ND	0.0125	12/17/2009
PFDS- Perfluorodecanesulfonate	ND	0.0125	12/16/2009
MeFOSAA- 2(N-Methylperfluorooctanesulfonamido) acetic acid	ND	0.0125	12/16/2009
EtFOSAA- 2(N-Ethylperfluorooctanesulfonamido) acetic acid	ND	0.0125	12/16/2009
MeFOSE- N-Methylperfluorooctane sulfonamidoethanol	ND	0.0125	12/16/2009
EtFOSE- N-Ethylperfluorooctane sulfonamidoethanol	ND	0.0125	12/16/2009

ND = Not detected = Response is below the LOD of 0.0125 ng/mL (ppb).

NQ = Not quantifiable = Response is between the LOD and the LOQ of 0.0250 ng/mL (ppb).



3058 Research Drive
State College, Pennsylvania 16801 USA
Telephone: 814.272.1039
Fax: 814.272.1019

Analytical Report

Fluorotelomer Analysis by GC/MS

Sample ID: L0019716-6; 17-3

Date of Extraction: 1/7/2010

Date Analyzed: 1/7/2010

Analyte	Amount Found (ng/mL)	LOQ (ng/mL)
7-2s FTOH	ND	1.00
6-2 FTOH	ND	1.00
8-2 FTOH	ND	1.00
10-2 FTOH	ND	1.00



Fluorochemical Residues in Water Samples By LC/MS/MS

Sample ID: L0019716-7; 18-1
Date of Extraction: 12/16/2009

Analyte	Amount Found (ng/mL)	LOD (ng/mL)	Analysis Date
C4 Acid- Perfluorobutyric Acid	ND	0.0125	12/16/2009
C5 Acid- Perfluoropentanoic Acid	ND	0.0125	12/16/2009
C6 Acid- Perfluorohexanoic Acid	ND	0.0125	12/16/2009
C7 Acid- Perfluoroheptanoic Acid	ND	0.0125	12/16/2009
C8 Acid- Perfluorooctanoic Acid	ND	0.0125	12/16/2009
C9 Acid- Perfluorononanoic Acid	ND	0.0125	12/16/2009
C10 Acid- Perfluorodecanoic Acid	ND	0.0125	12/16/2009
C11 Acid- Perfluoroundecanoic Acid	ND	0.0125	12/16/2009
C12 Acid- Perfluorododecanoic Acid	ND	0.0125	12/16/2009
C13 Acid- Perfluorotridecanoic Acid	ND	0.0125	12/16/2009
C14 Acid- Perfluorotetradecanoic Acid	ND	0.0125	12/16/2009
PFBS- Perfluorobutanesulfonate	NQ	0.0125	12/17/2009
PFHS- Perfluorohexanesulfonate	0.0350	0.0125	12/17/2009
PFOS- Perfluorooctanesulfonate	7.43	0.0125	12/17/2009
FOSA- Perfluorooctane sulfonamide	ND	0.0125	12/17/2009
PFHpS- Perfluoroheptanesulfonate	ND	0.0125	12/17/2009
PFDS- Perfluorodecanesulfonate	ND	0.0125	12/16/2009
MeFOSAA- 2(N-Methylperfluorooctanesulfonamido) acetic acid	ND	0.0125	12/16/2009
EtFOSAA- 2(N-Ethylperfluorooctanesulfonamido) acetic acid	ND	0.0125	12/16/2009
MeFOSE- N-Methylperfluorooctane sulfonamidoethanol	ND	0.0125	12/16/2009
EtFOSE- N-Ethylperfluorooctane sulfonamidoethanol	ND	0.0125	12/16/2009

ND = Not detected = Response is below the LOD of 0.0125 ng/mL (ppb).

NQ = Not quantifiable = Response is between the LOD and the LOQ of 0.0250 ng/mL (ppb).



3058 Research Drive
State College, Pennsylvania 16801 USA
Telephone: 814.272.1039
Fax: 814.272.1019

Analytical Report

Fluorotelomer Analysis by GC/MS

Sample ID: L0019716-9; 18-3

Date of Extraction: 1/12/2010

Date Analyzed: 1/12/2010

Analyte	Amount Found (ng/mL)	LOQ (ng/mL)
7-2s FTOH	ND	1.00
6-2 FTOH	ND	1.00
8-2 FTOH	ND	1.00
10-2 FTOH	ND	1.00

Fluorochemical Residues in Water Samples By LC/MS/MS

Sample ID: L0019716-10; 19-1
Date of Extraction: 12/16/2009

Analyte	Amount Found (ng/mL)	LOD (ng/mL)	Analysis Date
C4 Acid- Perfluorobutyric Acid	NQ	0.0125	12/16/2009
C5 Acid- Perfluoropentanoic Acid	0.0305	0.0125	12/16/2009
C6 Acid- Perfluorohexanoic Acid	0.207	0.0125	12/16/2009
C7 Acid- Perfluoroheptanoic Acid	0.0267	0.0125	12/16/2009
C8 Acid- Perfluorooctanoic Acid	0.101	0.0125	12/16/2009
C9 Acid- Perfluorononanoic Acid	0.144	0.0125	12/16/2009
C10 Acid- Perfluorodecanoic Acid	0.0452	0.0125	12/16/2009
C11 Acid- Perfluoroundecanoic Acid	ND	0.0125	12/16/2009
C12 Acid- Perfluorododecanoic Acid	ND	0.0125	12/16/2009
C13 Acid- Perfluorotridecanoic Acid	ND	0.0125	12/16/2009
C14 Acid- Perfluorotetradecanoic Acid	ND	0.0125	12/16/2009
PFBS- Perfluorobutanesulfonate	ND	0.0125	12/17/2009
PFHS- Perfluorohexanesulfonate	ND	0.0125	12/17/2009
PFOS- Perfluorooctanesulfonate	0.0404	0.0125	12/17/2009
FOSA- Perfluorooctane sulfonamide	ND	0.0125	12/17/2009
PFHpS- Perfluoroheptanesulfonate	ND	0.0125	12/17/2009
PFDS- Perfluorodecanesulfonate	ND	0.0125	12/16/2009
MeFOSAA- 2(N-Methylperfluorooctanesulfonamido) acetic acid	ND	0.0125	12/16/2009
EtFOSAA- 2(N-Ethylperfluorooctanesulfonamido) acetic acid	NQ	0.0125	12/16/2009
MeFOSE- N-Methylperfluorooctane sulfonamidoethanol	0.0462	0.0125	12/16/2009
EtFOSE- N-Ethylperfluorooctane sulfonamidoethanol	ND	0.0125	12/16/2009

ND = Not detected = Response is below the LOD of 0.0125 ng/mL (ppb).

NQ = Not quantifiable = Response is between the LOD and the LOQ of 0.0250 ng/mL (ppb).



3058 Research Drive
State College, Pennsylvania 16801 USA
Telephone: 814.272.1039
Fax: 814.272.1019

Analytical Report

Fluorotelomer Analysis by GC/MS

Sample ID: L0019716-12; 19-3

Date of Extraction: 1/12/2010

Date Analyzed: 1/12/2010

Analyte	Amount Found (ng/mL)	LOQ (ng/mL)
7-2s FTOH	ND	1.00
6-2 FTOH	ND	1.00
8-2 FTOH	ND	1.00
10-2 FTOH	ND	1.00



Fluorochemical Residues in Water Samples By LC/MS/MS

Sample ID: L0019716-13; 20-1

Date of Extraction: 12/16/2009

Analyte	Amount Found (ng/mL)	LOD (ng/mL)	Analysis Date
C4 Acid- Perfluorobutyric Acid	NQ	0.0125	12/16/2009
C5 Acid- Perfluoropentanoic Acid	0.0784	0.0125	12/16/2009
C6 Acid- Perfluorohexanoic Acid	0.159	0.0125	12/16/2009
C7 Acid- Perfluoroheptanoic Acid	0.0938	0.0125	12/16/2009
C8 Acid- Perfluorooctanoic Acid	0.277	0.0125	12/16/2009
C9 Acid- Perfluorononanoic Acid	0.0451	0.0125	12/16/2009
C10 Acid- Perfluorodecanoic Acid	0.0594	0.0125	12/16/2009
C11 Acid- Perfluoroundecanoic Acid	ND	0.0125	12/16/2009
C12 Acid- Perfluorododecanoic Acid	ND	0.0125	12/16/2009
C13 Acid- Perfluorotridecanoic Acid	ND	0.0125	12/16/2009
C14 Acid- Perfluorotetradecanoic Acid	ND	0.0125	12/16/2009
PFBS- Perfluorobutanesulfonate	ND	0.0125	12/17/2009
PFHS- Perfluorohexanesulfonate	ND	0.0125	12/17/2009
PFOS- Perfluorooctanesulfonate	0.107	0.0125	12/17/2009
FOSA- Perfluorooctane sulfonamide	ND	0.0125	12/17/2009
PFHpS- Perfluoroheptanesulfonate	ND	0.0125	12/17/2009
PFDS- Perfluorodecanesulfonate	ND	0.0125	12/16/2009
MeFOSAA- 2(N-Methylperfluorooctanesulfonamido) acetic acid	0.146	0.0125	12/16/2009
EtFOSAA- 2(N-Ethylperfluorooctanesulfonamido) acetic acid	ND	0.0125	12/16/2009
MeFOSE- N-Methylperfluorooctane sulfonamidoethanol	0.0555	0.0125	12/16/2009
EtFOSE- N-Ethylperfluorooctane sulfonamidoethanol	ND	0.0125	12/16/2009

ND = Not detected = Response is below the LOD of 0.0125 ng/mL (ppb).

NQ = Not quantifiable = Response is between the LOD and the LOQ of 0.0250 ng/mL (ppb).

Fluorotelomer Analysis by GC/MS**Sample ID: L0019716-15; 20-3****Date of Extraction: 1/7/2010****Date Analyzed: 1/7/2010**

Analyte	Amount Found (ng/mL)	LOQ (ng/mL)
7-2s FTOH	ND	1.00
6-2 FTOH	1.40	1.00
8-2 FTOH	ND	1.00
10-2 FTOH	ND	1.00



Fluorochemical Residues in Water Samples By LC/MS/MS

Sample ID: L0019716-16; 21-1
Date of Extraction: 12/16/2009

Analyte	Amount Found (ng/mL)	LOD (ng/mL)	Analysis Date
C4 Acid- Perfluorobutyric Acid	NQ	0.0125	12/16/2009
C5 Acid- Perfluoropentanoic Acid	0.0394	0.0125	12/16/2009
C6 Acid- Perfluorohexanoic Acid	0.0696	0.0125	12/16/2009
C7 Acid- Perfluoroheptanoic Acid	NQ	0.0125	12/16/2009
C8 Acid- Perfluorooctanoic Acid	0.0389	0.0125	12/16/2009
C9 Acid- Perfluorononanoic Acid	ND	0.0125	12/16/2009
C10 Acid- Perfluorodecanoic Acid	ND	0.0125	12/16/2009
C11 Acid- Perfluoroundecanoic Acid	ND	0.0125	12/16/2009
C12 Acid- Perfluorododecanoic Acid	ND	0.0125	12/16/2009
C13 Acid- Perfluorotridecanoic Acid	ND	0.0125	12/16/2009
C14 Acid- Perfluorotetradecanoic Acid	ND	0.0125	12/16/2009
PFBS- Perfluorobutanesulfonate	NQ	0.0125	12/17/2009
PFHS- Perfluorohexanesulfonate	ND	0.0125	12/17/2009
PFOS- Perfluorooctanesulfonate	0.233	0.0125	12/17/2009
FOSA- Perfluorooctane sulfonamide	ND	0.0125	12/17/2009
PFHpS- Perfluoroheptanesulfonate	ND	0.0125	12/17/2009
PFDS- Perfluorodecanesulfonate	ND	0.0125	12/16/2009
MeFOSAA- 2(N-Methylperfluorooctanesulfonamido) acetic acid	ND	0.0125	12/16/2009
EtFOSAA- 2(N-Ethylperfluorooctanesulfonamido) acetic acid	ND	0.0125	12/16/2009
MeFOSE- N-Methylperfluorooctane sulfonamidoethanol	ND	0.0125	12/16/2009
EtFOSE- N-Ethylperfluorooctane sulfonamidoethanol	ND	0.0125	12/16/2009

ND = Not detected = Response is below the LOD of 0.0125 ng/mL (ppb).

NQ = Not quantifiable = Response is between the LOD and the LOQ of 0.0250 ng/mL (ppb).



3058 Research Drive
State College, Pennsylvania 16801 USA
Telephone: 814.272.1039
Fax: 814.272.1019

Analytical Report

Fluorotelomer Analysis by GC/MS

Sample ID: L0019716-18; 21-3

Date of Extraction: 1/12/2010

Date Analyzed: 1/12/2010

Analyte	Amount Found (ng/mL)	LOQ (ng/mL)
7-2s FTOH	ND	1.00
6-2 FTOH	2.60	1.00
8-2 FTOH	ND	1.00
10-2 FTOH	ND	1.00



Fluorochemical Residues in Water Samples By LC/MS/MS

Sample ID: L0019716-19; 22-1

Date of Extraction: 12/16/2009

Analyte	Amount Found (ng/mL)	LOD (ng/mL)	Analysis Date
C4 Acid- Perfluorobutyric Acid	NQ	0.0125	12/16/2009
C5 Acid- Perfluoropentanoic Acid	NQ	0.0125	12/16/2009
C6 Acid- Perfluorohexanoic Acid	ND	0.0125	12/16/2009
C7 Acid- Perfluoroheptanoic Acid	ND	0.0125	12/16/2009
C8 Acid- Perfluorooctanoic Acid	0.137	0.0125	12/16/2009
C9 Acid- Perfluorononanoic Acid	ND	0.0125	12/16/2009
C10 Acid- Perfluorodecanoic Acid	ND	0.0125	12/16/2009
C11 Acid- Perfluoroundecanoic Acid	ND	0.0125	12/16/2009
C12 Acid- Perfluorododecanoic Acid	ND	0.0125	12/16/2009
C13 Acid- Perfluorotridecanoic Acid	ND	0.0125	12/16/2009
C14 Acid- Perfluorotetradecanoic Acid	ND	0.0125	12/16/2009
PFBS- Perfluorobutanesulfonate	0.0347	0.0125	12/17/2009
PFHS- Perfluorohexanesulfonate	ND	0.0125	12/17/2009
PFOS- Perfluorooctanesulfonate	1.90	0.0125	12/17/2009
FOSA- Perfluorooctane sulfonamide	ND	0.0125	12/17/2009
PFHpS- Perfluoroheptanesulfonate	ND	0.0125	12/17/2009
PFDS- Perfluorodecanesulfonate	ND	0.0125	12/16/2009
MeFOSAA- 2(N-Methylperfluorooctanesulfonamido) acetic acid	ND	0.0125	12/16/2009
EtFOSAA- 2(N-Ethylperfluorooctanesulfonamido) acetic acid	ND	0.0125	12/16/2009
MeFOSE- N-Methylperfluorooctane sulfonamidoethanol	0.348	0.0125	12/16/2009
EtFOSE- N-Ethylperfluorooctane sulfonamidoethanol	ND	0.0125	12/16/2009

ND = Not detected = Response is below the LOD of 0.0125 ng/mL (ppb).

NQ = Not quantifiable = Response is between the LOD and the LOQ of 0.0250 ng/mL (ppb).



3058 Research Drive
State College, Pennsylvania 16801 USA
Telephone: 814.272.1039
Fax: 814.272.1019

Analytical Report

Fluorotelomer Analysis by GC/MS

Sample ID: L0019716-21; 22-3

Date of Extraction: 1/7/2010

Date Analyzed: 1/7/2010

Analyte	Amount Found (ng/mL)	LOQ (ng/mL)
7-2s FTOH	ND	1.00
6-2 FTOH	ND	1.00
8-2 FTOH	ND	1.00
10-2 FTOH	ND	1.00



Fluorochemical Residues in Water Samples By LC/MS/MS

Sample ID: L0019716-22; 23-1
Date of Extraction: 12/16/2009

Analyte	Amount Found (ng/mL)	LOD (ng/mL)	Analysis Date
C4 Acid- Perfluorobutyric Acid	ND	0.0125	12/16/2009
C5 Acid- Perfluoropentanoic Acid	NQ	0.0125	12/16/2009
C6 Acid- Perfluorohexanoic Acid	0.0296	0.0125	12/16/2009
C7 Acid- Perfluoroheptanoic Acid	0.0470	0.0125	12/16/2009
C8 Acid- Perfluorooctanoic Acid	0.0900	0.0125	12/16/2009
C9 Acid- Perfluorononanoic Acid	0.0368	0.0125	12/16/2009
C10 Acid- Perfluorodecanoic Acid	ND	0.0125	12/16/2009
C11 Acid- Perfluoroundecanoic Acid	NQ	0.0125	12/16/2009
C12 Acid- Perfluorododecanoic Acid	ND	0.0125	12/16/2009
C13 Acid- Perfluorotridecanoic Acid	ND	0.0125	12/16/2009
C14 Acid- Perfluorotetradecanoic Acid	ND	0.0125	12/16/2009
PFBS- Perfluorobutanesulfonate	ND	0.0125	12/17/2009
PFHS- Perfluorohexanesulfonate	ND	0.0125	12/17/2009
PFOS- Perfluorooctanesulfonate	1.81	0.0125	12/17/2009
FOSA- Perfluorooctane sulfonamide	ND	0.0125	12/17/2009
PFHpS- Perfluoroheptanesulfonate	ND	0.0125	12/17/2009
PFDS- Perfluorodecanesulfonate	ND	0.0125	12/16/2009
MeFOSAA- 2(N-Methylperfluorooctanesulfonamido) acetic acid	ND	0.0125	12/16/2009
EtFOSAA- 2(N-Ethylperfluorooctanesulfonamido) acetic acid	ND	0.0125	12/16/2009
MeFOSE- N-Methylperfluorooctane sulfonamidoethanol	ND	0.0125	12/16/2009
EtFOSE- N-Ethylperfluorooctane sulfonamidoethanol	ND	0.0125	12/16/2009

ND = Not detected = Response is below the LOD of 0.0125 ng/mL (ppb).

NQ = Not quantifiable = Response is between the LOD and the LOQ of 0.0250 ng/mL (ppb).



3058 Research Drive
State College, Pennsylvania 16801 USA
Telephone: 814.272.1039
Fax: 814.272.1019

Analytical Report

Fluorotelomer Analysis by GC/MS

Sample ID: L0019716-24; 23-3

Date of Extraction: 1/7/2010

Date Analyzed: 1/7/2010

Analyte	Amount Found (ng/mL)	LOQ (ng/mL)
7-2s FTOH	ND	1.00
6-2 FTOH	ND	1.00
8-2 FTOH	ND	1.00
10-2 FTOH	ND	1.00



MPI
RESEARCH

3058 Research Drive
State College, Pennsylvania 16801 USA
Telephone: 814.272.1039
Fax: 814.272.1019

Analytical Report

Recovery Summary of 8:1 FTOH in Water Samples

Client Sample ID	MPI Sample ID	Amount Spiked (ng/mL, ppb)	Amount Recovered (ng/mL, ppb)	Recovery (%)
N/A	Method Blank (1/7/10)	5.00	3.30	66
N/A	LCS (1/7/10)	5.00	3.00	60
N/A	LCSD (1/7/10)	5.00	3.50	70
N/A	Method Blank (1/12/10)	5.00	3.60	72
N/A	LCS (1/12/10)	5.00	3.30	66
N/A	LCSD (1/12/10)	5.00	3.70	74
N/A	Method Blank (1/19/10)	5.00	3.30	66
N/A	LCS (1/19/10)	5.00	2.90	58
N/A	LCSD (1/19/10)	5.00	3.40	68
16-3	L19716-3	5.00	2.30	46
17-3	L19716-6	5.00	3.50	70
18-3	L19716-9	5.00	1.00	20
19-3	L19716-12	5.00	1.10	22
20-3	L19716-15	5.00	2.70	54
21-3	L19716-18	5.00	1.30	26
22-3	L19716-21	5.00	3.40	68
23-3	L19716-24	5.00	3.40	68
24-3	L19716-27	5.00	2.80	56
2D3	L19716-30	5.00	1.00	20
TRIP (FTOH)	L19716-31	5.00	2.80	56

Fluorochemical Residues in Water Samples By LC/MS/MS

Sample ID: L0019716-25; 24-1

Date of Extraction: 12/16/2009

Analyte	Amount Found (ng/mL)	LOD (ng/mL)	Analysis Date
C4 Acid- Perfluorobutyric Acid	ND	0.0125	12/16/2009
C5 Acid- Perfluoropentanoic Acid	ND	0.0125	12/16/2009
C6 Acid- Perfluorohexanoic Acid	ND	0.0125	12/16/2009
C7 Acid- Perfluoroheptanoic Acid	ND	0.0125	12/16/2009
C8 Acid- Perfluorooctanoic Acid	0.0258	0.0125	12/16/2009
C9 Acid- Perfluorononanoic Acid	ND	0.0125	12/16/2009
C10 Acid- Perfluorodecanoic Acid	ND	0.0125	12/16/2009
C11 Acid- Perfluoroundecanoic Acid	ND	0.0125	12/16/2009
C12 Acid- Perfluorododecanoic Acid	ND	0.0125	12/16/2009
C13 Acid- Perfluorotridecanoic Acid	ND	0.0125	12/16/2009
C14 Acid- Perfluorotetradecanoic Acid	ND	0.0125	12/16/2009
PFBS- Perfluorobutanesulfonate	ND	0.0125	12/17/2009
PFHS- Perfluorohexanesulfonate	ND	0.0125	12/17/2009
PFOS- Perfluorooctanesulfonate	0.0618	0.0125	12/17/2009
FOSA- Perfluorooctane sulfonamide	ND	0.0125	12/17/2009
PFHpS- Perfluoroheptanesulfonate	ND	0.0125	12/17/2009
PFDS- Perfluorodecanesulfonate	ND	0.0125	12/16/2009
MeFOSAA- 2(N-Methylperfluorooctanesulfonamido) acetic acid	ND	0.0125	12/16/2009
EtFOSAA- 2(N-Ethylperfluorooctanesulfonamido) acetic acid	ND	0.0125	12/16/2009
MeFOSE- N-Methylperfluorooctane sulfonamidoethanol	ND	0.0125	12/16/2009
EtFOSE- N-Ethylperfluorooctane sulfonamidoethanol	ND	0.0125	12/16/2009

ND = Not detected = Response is below the LOD of 0.0125 ng/mL (ppb).

NQ = Not quantifiable = Response is between the LOD and the LOQ of 0.0250 ng/mL (ppb).



MPI
RESEARCH

3058 Research Drive
State College, Pennsylvania 16801 USA
Telephone: 814.272.1039
Fax: 814.272.1019

Analytical Report

Fluorotelomer Analysis by GC/MS

Sample ID: L0019716-27; 24-3

Date of Extraction: 1/7/2010

Date Analyzed: 1/7/2010

Analyte	Amount Found (ng/mL)	LOQ (ng/mL)
7-2s FTOH	ND	1.00
6-2 FTOH	ND	1.00
8-2 FTOH	ND	1.00
10-2 FTOH	ND	1.00

Fluorochemical Residues in Water Samples By LC/MS/MS

Sample ID: L0019716-28; 3D1
Date of Extraction: 12/16/2009

Analyte	Amount Found (ng/mL)	LOD (ng/mL)	Analysis Date
C4 Acid- Perfluorobutyric Acid	ND	0.0125	12/16/2009
C5 Acid- Perfluoropentanoic Acid	ND	0.0125	12/16/2009
C6 Acid- Perfluorohexanoic Acid	ND	0.0125	12/16/2009
C7 Acid- Perfluoroheptanoic Acid	ND	0.0125	12/16/2009
C8 Acid- Perfluorooctanoic Acid	ND	0.0125	12/16/2009
C9 Acid- Perfluorononanoic Acid	NQ	0.0125	12/16/2009
C10 Acid- Perfluorodecanoic Acid	ND	0.0125	12/16/2009
C11 Acid- Perfluoroundecanoic Acid	ND	0.0125	12/16/2009
C12 Acid- Perfluorododecanoic Acid	ND	0.0125	12/16/2009
C13 Acid- Perfluorotridecanoic Acid	ND	0.0125	12/16/2009
C14 Acid- Perfluorotetradecanoic Acid	ND	0.0125	12/16/2009
PFBS- Perfluorobutanesulfonate	NQ	0.0125	12/17/2009
PFHS- Perfluorohexanesulfonate	0.0323	0.0125	12/17/2009
PFOS- Perfluorooctanesulfonate	7.29	0.0125	12/17/2009
FOSA- Perfluorooctane sulfonamide	ND	0.0125	12/17/2009
PFHpS- Perfluoroheptanesulfonate	ND	0.0125	12/17/2009
PFDS- Perfluorodecanesulfonate	ND	0.0125	12/16/2009
MeFOSAA- 2(N-Methylperfluorooctanesulfonamido) acetic acid	ND	0.0125	12/16/2009
EtFOSAA- 2(N-Ethylperfluorooctanesulfonamido) acetic acid	ND	0.0125	12/16/2009
MeFOSE- N-Methylperfluorooctane sulfonamidoethanol	NQ	0.0125	12/16/2009
EtFOSE- N-Ethylperfluorooctane sulfonamidoethanol	ND	0.0125	12/16/2009

ND = Not detected = Response is below the LOD of 0.0125 ng/mL (ppb).

NQ = Not quantifiable = Response is between the LOD and the LOQ of 0.0250 ng/mL (ppb).



3058 Research Drive
State College, Pennsylvania 16801 USA
Telephone: 814.272.1039
Fax: 814.272.1019

Analytical Report

Fluorotelomer Analysis by GC/MS

Sample ID: L0019716-30; 2D3

Date of Extraction: 1/12/2010

Date Analyzed: 1/12/2010

Analyte	Amount Found (ng/mL)	LOQ (ng/mL)
7-2s FTOH	ND	1.00
6-2 FTOH	ND	1.00
8-2 FTOH	ND	1.00
10-2 FTOH	ND	1.00

Fluorochemical Residues in Water Samples By LC/MS/MS

Sample ID: L0019716-32; TRIP (PFCs)

Date of Extraction: 12/16/2009

Analyte	Amount Found (ng/mL)	LOD (ng/mL)	Analysis Date
C4 Acid- Perfluorobutyric Acid	ND	0.0125	12/16/2009
C5 Acid- Perfluoropentanoic Acid	ND	0.0125	12/16/2009
C6 Acid- Perfluorohexanoic Acid	ND	0.0125	12/16/2009
C7 Acid- Perfluoroheptanoic Acid	ND	0.0125	12/16/2009
C8 Acid- Perfluorooctanoic Acid	ND	0.0125	12/16/2009
C9 Acid- Perfluorononanoic Acid	ND	0.0125	12/16/2009
C10 Acid- Perfluorodecanoic Acid	ND	0.0125	12/16/2009
C11 Acid- Perfluoroundecanoic Acid	ND	0.0125	12/16/2009
C12 Acid- Perfluorododecanoic Acid	ND	0.0125	12/16/2009
C13 Acid- Perfluorotridecanoic Acid	ND	0.0125	12/16/2009
C14 Acid- Perfluorotetradecanoic Acid	ND	0.0125	12/16/2009
PFBS- Perfluorobutanesulfonate	ND	0.0125	12/17/2009
PFHS- Perfluorohexanesulfonate	ND	0.0125	12/17/2009
PFOS- Perfluorooctanesulfonate	ND	0.0125	12/17/2009
FOSA- Perfluorooctane sulfonamide	ND	0.0125	12/17/2009
PFHpS- Perfluoroheptanesulfonate	ND	0.0125	12/17/2009
PFDS- Perfluorodecanesulfonate	ND	0.0125	12/16/2009
MeFOSAA- 2(N-Methylperfluorooctanesulfonamido) acetic acid	ND	0.0125	12/16/2009
EtFOSAA- 2(N-Ethylperfluorooctanesulfonamido) acetic acid	ND	0.0125	12/16/2009
MeFOSE- N-Methylperfluorooctane sulfonamidoethanol	ND	0.0125	12/16/2009
EtFOSE- N-Ethylperfluorooctane sulfonamidoethanol	ND	0.0125	12/16/2009

ND = Not detected = Response is below the LOD of 0.0125 ng/mL (ppb).

NQ = Not quantifiable = Response is between the LOD and the LOQ of 0.0250 ng/mL (ppb).



MPI
RESEARCH

3058 Research Drive
State College, Pennsylvania 16801 USA
Telephone: 814.272.1039
Fax: 814.272.1019

Analytical Report

Fluorotelomer Analysis by GC/MS

Sample ID: L0019716-31; TRIP (FTOH)

Date of Extraction: 1/19/2010

Date Analyzed: 1/19/2010

Analyte	Amount Found (ng/mL)	LOQ (ng/mL)
7-2s FTOH	ND	1.00
6-2 FTOH	ND	1.00
8-2 FTOH	ND	1.00
10-2 FTOH	ND	1.00

Recovery Summary of Fluorochemical Residues In Water Samples

Sample Description	Amount Spiked (ng/mL)	C4 Acid			C5 Acid			C6 Acid		
		Amt Found In Sample (ng/mL)	Amount Recovered (ng/mL)	Recovery (%)	Amt Found In Sample (ng/mL)	Amount Recovered (ng/mL)	Recovery (%)	Amt Found In Sample (ng/mL)	Amount Recovered (ng/mL)	Recovery (%)
Reagent Spike A 0.1 ng/mL	0.100	ND	0.101	101	ND	0.104	104	ND	0.100	100
Reagent Spike B 1.0 ng/mL	1.00	ND	1.09	109	ND	1.10	110	ND	1.20	120
16-1 Spike C (L19716-1 Spk C, 1.0 ng/mL Lab Spike)	1.00	NQ	0.996	100	0.102	1.21	111	0.238	1.33	109

Sample Description	Amount Spiked (ng/mL)	C7 Acid			C8 Acid			C9 Acid		
		Amt Found In Sample (ng/mL)	Amount Recovered (ng/mL)	Recovery (%)	Amt Found In Sample (ng/mL)	Amount Recovered (ng/mL)	Recovery (%)	Amt Found In Sample (ng/mL)	Amount Recovered (ng/mL)	Recovery (%)
Reagent Spike A 0.1 ng/mL	0.100	ND	0.103	103	ND	0.100	100	ND	0.108	108
Reagent Spike B 1.0 ng/mL	1.00	ND	1.18	118	ND	1.11	111	ND	1.17	117
16-1 Spike C (L19716-1 Spk C, 1.0 ng/mL Lab Spike)	1.00	0.0289	1.07	104	0.0320	1.13	110	ND	1.14	114

Sample Description	Amount Spiked (ng/mL)	C10 Acid			C11 Acid			C12 Acid		
		Amt Found In Sample (ng/mL)	Amount Recovered (ng/mL)	Recovery (%)	Amt Found In Sample (ng/mL)	Amount Recovered (ng/mL)	Recovery (%)	Amt Found In Sample (ng/mL)	Amount Recovered (ng/mL)	Recovery (%)
Reagent Spike A 0.1 ng/mL	0.100	ND	0.106	106	ND	0.124	124	ND	0.107	107
Reagent Spike B 1.0 ng/mL	1.00	ND	1.12	112	ND	1.25	125	ND	1.21	121
16-1 Spike C (L19716-1 Spk C, 1.0 ng/mL Lab Spike)	1.00	0.0255	1.16	113	ND	1.33	133^	ND	1.14	114

Sample Description	Amount Spiked (ng/mL)	C13 Acid			C14 Acid		
		Amt Found In Sample (ng/mL)	Amount Recovered (ng/mL)	Recovery (%)	Amt Found In Sample (ng/mL)	Amount Recovered (ng/mL)	Recovery (%)
Reagent Spike A 0.1 ng/mL	0.100	ND	0.112	112	ND	0.119	119
Reagent Spike B 1.0 ng/mL	1.00	ND	1.24	124	ND	1.21	121
16-1 Spike C (L19716-1 Spk C, 1.0 ng/mL Lab Spike)	1.00	ND	1.59	159^	ND	1.41	141^

ND = Not detected = Response is below the LOD of 0.0125 ng/mL.

NQ = Not quantifiable = Response is between the LOD and the LOQ of 0.0250 ng/mL.

^ Laboratory Matrix Spike recovery is outside the acceptance criteria of 70 - 130%. The recovery is bias high and the result is either non detect or non quantifiable, therefore, the data is considered acceptable.

Recovery Summary of Fluorochemical Residues in Water Samples

Sample Description	Amount Spiked (ng/mL)	PFBS			PFHS			PFOS		
		Amt Found In Sample (ng/mL)	Amount Recovered (ng/mL)	Recovery (%)	Amt Found In Sample (ng/mL)	Amount Recovered (ng/mL)	Recovery (%)	Amt Found In Sample (ng/mL)	Amount Recovered (ng/mL)	Recovery (%)
Reagent Spike A 0.1 ng/mL	0.100	ND	0.0947	95	ND	0.104	104	ND	0.104	104
Reagent Spike B 1.0 ng/mL	1.00	ND	1.08	108	ND	1.25	125	ND	1.18	118
16-1 Spike C (L19716-1 Spk C, 1.0 ng/mL Lab Spike)	1.00	ND	0.990	99	ND	1.05	105	0.162	1.30	114

Sample Description	Amount Spiked (ng/mL)	FOSA			PFHpS			PFDS		
		Amt Found In Sample (ng/mL)	Amount Recovered (ng/mL)	Recovery (%)	Amt Found In Sample (ng/mL)	Amount Recovered (ng/mL)	Recovery (%)	Amt Found In Sample (ng/mL)	Amount Recovered (ng/mL)	Recovery (%)
Reagent Spike A 0.1 ng/mL	0.100	ND	0.101	101	ND	0.107	107	ND	0.103	103
Reagent Spike B 1.0 ng/mL	1.00	ND	1.15	115	ND	1.13	113	ND	1.14	114
16-1 Spike C (L19716-1 Spk C, 1.0 ng/mL Lab Spike)	1.00	ND	1.07	107	ND	1.02	102	ND	1.05	105

ND = Not detected = Response is below the LOD of 0.0125 ng/mL.

NQ = Not quantifiable = Response is between the LOD and the LOQ of 0.0250 ng/mL.

Recovery Summary of Fluorochemical Residues in Water Samples

Sample Description	MeFOSAA				EtFOSAA		
	Amount Spiked (ng/mL)	Amt Found in Sample (ng/mL)	Amount Recovered (ng/mL)	Recovery (%)	Amt Found in Sample (ng/mL)	Amount Recovered (ng/mL)	Recovery (%)
Reagent Spike A 0.1 ng/mL	0.100	ND	0.100	100	ND	0.103	103
Reagent Spike B 1.0 ng/mL	1.00	ND	1.02	102	ND	1.01	101
16-1 Spike C (L19716-1 Spk C, 0.1 ng/mL Lab Spike)	0.100	ND	0.113	113	ND	0.129	129

Sample Description	MeFOSE				EtFOSE		
	Amount Spiked (ng/mL)	Amt Found in Sample (ng/mL)	Amount Recovered (ng/mL)	Recovery (%)	Amt Found in Sample (ng/mL)	Amount Recovered (ng/mL)	Recovery (%)
Reagent Spike A 0.1 ng/mL	0.100	ND	0.106	106	ND	0.112	112
Reagent Spike B 1.0 ng/mL	1.00	ND	0.989	99	ND	1.01	101
16-1 Spike C (L19716-1 Spk C, 0.1 ng/mL Lab Spike)	0.100	ND	0.110	110	ND	0.0829	83

ND = Not detected = Response is below the LOD of 0.0125 ng/mL.

NQ = Not quantifiable = Response is between the LOD and the LOQ of 0.0250 ng/mL.



MPI
RESEARCH

3058 Research Drive
State College, Pennsylvania 16801 USA
Telephone: 814.272.1039
Fax: 814.272.1019

Analytical Report

Recovery Summary of ¹³C PFOA (m+4) in Water Samples

Client Sample ID	MPI Sample ID	Amount Spiked (ng/mL, ppb)	Amount Recovered (ng/mL, ppb)	Recovery (%)
N/A	Reagent Control	1.00	1.04	104
N/A	Reagent Spike A	0.100	0.104	104
N/A	Reagent Spike B	1.00	1.19	119
16-1 Spike C	L19716-1 Spk C	1.00	1.05	105
16-1	L19716-1	1.00	0.979	98
17-1	L19716-4	1.00	0.991	99
18-1	L19716-7	1.00	1.08	108
19-1	L19716-10	1.00	0.984	98
20-1	L19716-13	1.00	1.08	108
21-1	L19716-16	1.00	1.08	108
22-1	L19716-19	1.00	0.983	98
23-1	L19716-22	1.00	1.08	108
24-1	L19716-25	1.00	1.08	108
3D1	L19716-28	1.00	0.983	98
TRIP (PFCs)	L19716-32	1.00	1.01	101

Fluorotelomer Analysis by GC/MS

Sample ID: Method Blank
Date of Extraction: 1/7/2010
Date Analyzed: 1/7/2010

Analyte	Amount Found (ng/mL)	LOQ (ng/mL)
7-2s FTOH	ND	1.00
6-2 FTOH	ND	1.00
8-2 FTOH	ND	1.00
10-2 FTOH	ND	1.00



MPI
RESEARCH

3058 Research Drive
State College, Pennsylvania 16801 USA
Telephone: 814.272.1039
Fax: 814.272.1019

Analytical Report

Fluorotelomer Analysis by GC/MS

Sample ID: Lab Contol Spike

Date of Extraction: 1/7/2010

Date Analyzed: 1/7/2010

Analyte	Amount Found (ng/mL)	Amount Spiked (ng/mL)	LOQ (ng/mL)	% Recovery
7-2s FTOH	3.30	5.00	1.00	66
6-2 FTOH	7.90	9.88	1.00	80
8-2 FTOH	3.70	5.11	1.00	72
10-2 FTOH	3.60	5.03	1.00	72



3058 Research Drive
State College, Pennsylvania 16801 USA
Telephone: 814.272.1039
Fax: 814.272.1019

Analytical Report

Fluorotelomer Analysis by GC/MS

Sample ID: Lab Control Spike Duplicate

Date of Extraction: 1/7/2010

Date Analyzed: 1/7/2010

Analyte	Amount Found (ng/mL)	Amount Spiked (ng/mL)	LOQ (ng/mL)	% Recovery	% RPD
7-2s FTOH	3.70	5.00	1.00	74	11.4
6-2 FTOH	8.80	9.88	1.00	89	10.8
8-2 FTOH	4.20	5.11	1.00	82	12.7
10-2 FTOH	4.10	5.03	1.00	82	13.0



MPI
RESEARCH

3058 Research Drive
State College, Pennsylvania 16801 USA
Telephone: 814.272.1039
Fax: 814.272.1019

Analytical Report

Fluorotelomer Analysis by GC/MS

Sample ID: Method Blank
Date of Extraction: 1/12/2010
Date Analyzed: 1/12/2010

Analyte	Amount Found (ng/mL)	LOQ (ng/mL)
7-2s FTOH	ND	1.00
6-2 FTOH	ND	1.00
8-2 FTOH	ND	1.00
10-2 FTOH	ND	1.00



3058 Research Drive
State College, Pennsylvania 16801 USA
Telephone: 814.272.1039
Fax: 814.272.1019

Analytical Report

Fluorotelomer Analysis by GC/MS

Sample ID: Lab Control Spike

Date of Extraction: 1/12/2010

Date Analyzed: 1/13/2010

Analyte	Amount Found (ng/mL)	Amount Spiked (ng/mL)	LOQ (ng/mL)	% Recovery
7-2s FTOH	3.30	5.00	1.00	66
6-2 FTOH	8.20	9.88	1.00	83
8-2 FTOH	3.60	5.11	1.00	71
10-2 FTOH	3.60	5.03	1.00	72



MPI
RESEARCH

3058 Research Drive
State College, Pennsylvania 16801 USA
Telephone: 814.272.1039
Fax: 814.272.1019

Analytical Report

Fluorotelomer Analysis by GC/MS

Sample ID: Lab Control Spike Duplicate

Date of Extraction: 1/12/2010

Date Analyzed: 1/12/2010

Analyte	Amount Found (ng/mL)	Amount Spiked (ng/mL)	LOQ (ng/mL)	% Recovery	% RPD
7-2s FTOH	3.40	5.00	1.00	68	3.0
6-2 FTOH	8.50	9.88	1.00	86	3.6
8-2 FTOH	4.00	5.11	1.00	78	10.5
10-2 FTOH	4.00	5.03	1.00	80	10.5

Fluorotelomer Analysis by GC/MS

Sample ID: Method Blank
Date of Extraction: 1/19/2010
Date Analyzed: 1/19/2010

Analyte	Amount Found (ng/mL)	LOQ (ng/mL)
7-2s FTOH	ND	1.00
6-2 FTOH	ND	1.00
8-2 FTOH	ND	1.00
10-2 FTOH	ND	1.00



MPI
RESEARCH

3058 Research Drive
State College, Pennsylvania 16801 USA
Telephone: 814.272.1039
Fax: 814.272.1019

Analytical Report

Fluorotelomer Analysis by GC/MS

Sample ID: Lab Contol Spike

Date of Extraction: 1/19/2010

Date Analyzed: 1/19/2010

Analyte	Amount Found (ng/mL)	Amount Spiked (ng/mL)	LOQ (ng/mL)	% Recovery
7-2s FTOH	3.40	5.00	1.00	68
6-2 FTOH	3.40	4.94	1.00	69
8-2 FTOH	3.30	5.11	1.00	65
10-2 FTOH	3.40	5.03	1.00	68



3058 Research Drive
State College, Pennsylvania 16801 USA
Telephone: 814.272.1039
Fax: 814.272.1019

Analytical Report

Fluorotelomer Analysis by GC/MS

Sample ID: Lab Control Spike Duplicate

Date of Extraction: 1/19/2010

Date Analyzed: 1/19/2010

Analyte	Amount Found (ng/mL)	Amount Spiked (ng/mL)	LOQ (ng/mL)	% Recovery	% RPD
7-2s FTOH	3.70	5.00	1.00	74	8.8
6-2 FTOH	3.50	4.94	1.00	71	2.9
8-2 FTOH	3.30	5.11	1.00	65	0.0
10-2 FTOH	3.60	5.03	1.00	72	5.1