



**MPI**  
RESEARCH

3058 Research Drive  
State College, Pennsylvania 16801 USA  
Telephone: 814.272.1039  
Fax: 814.272.1019

---

## **Analytical Report**

**Fluorochemical Characterization of Aqueous Samples**

**MPI Report No. L0019874**

---

### ***Testing Laboratory***

MPI Research, Inc.  
3058 Research Drive  
State College, PA 16801

---

### ***Requester***

Jason Governo  
University of Georgia  
Driftmier Engineering Center  
Athens, GA 30602  
Phone: 678-794-6664

## 1 Introduction

Results are reported for the analysis of twelve water samples received at MPI Research from the University of Georgia. The MPI Research study number assigned to the project is L0019874. Table I lists the target analyte quantitated for the sample.

Table I. Target Analyte for Quantitation

| Compound Name                          | Acronym               | Analysis |
|--|-----------------------|----------|
| Perfluorobutyric Acid                  | C4 Acid               | LC/MS/MS |
| Perfluoropentanoic Acid                | C5 Acid               | LC/MS/MS |
| Perfluorohexanoic Acid                 | C6 Acid               | LC/MS/MS |
| Perfluoroheptanoic Acid                | C7 Acid               | LC/MS/MS |
| Perfluorooctanoic Acid                 | C8 Acid               | LC/MS/MS |
| Perfluorononanoic Acid                 | C9 Acid               | LC/MS/MS |
| Perfluorodecanoic Acid                 | C10 Acid              | LC/MS/MS |
| Perfluoroundecanoic Acid               | C11 Acid              | LC/MS/MS |
| Perfluorododecanoic Acid               | C12 Acid              | LC/MS/MS |
| Perfluorotridecanoic Acid              | C13 Acid              | LC/MS/MS |
| Perfluorotetradecanoic Acid            | C14 Acid              | LC/MS/MS |
| Perfluorobutanesulfonate               | C4 Sulfonate or PFBS  | LC/MS/MS |
| Perfluorohexanesulfonate               | C6 Sulfonate or PFHS  | LC/MS/MS |
| Perfluoroheptanesulfonate              | C7 Sulfonate or PFOS  | LC/MS/MS |
| Perfluorooctanesulfonate               | C8 Sulfonate or PFOS  | LC/MS/MS |
| Perfluorodecanesulfonate               | C10 Sulfonate or PFOS | LC/MS/MS |
| Perfluorooctanesulfonamide             | FOSA                  | LC/MS/MS |
| 2-(N-methylperfluorooctanesulfonamido) | MeFOSAA               | LC/MS/MS |
| 2-(N-ethylperfluorooctanesulfonamido)  | EtFOSAA               | LC/MS/MS |
| N-methylperfluorooctane                | MeFOSE                | LC/MS/MS |
| N-ethylperfluorooctane                 | EtFOSE                | LC/MS/MS |
| 6:2 Fluorotelomer alcohol              | 6:2 FTOH              | GC/MS    |
| 7:2 sFluorotelomer alcohol             | 7:2s FTOH             | GC/MS    |
| 8:2 Fluorotelomer alcohol              | 8:2 FTOH              | GC/MS    |
| 10:2 Fluorotelomer alcohol             | 10:2 FTOH             | GC/MS    |

## 2 Sample Receipt

Twelve water samples were received from the client cooled with wet ice on 12/23/09 and given the MPI Research login number of L0019874. The samples were stored refrigerated from receipt until analysis. Chain-of-custody information is presented in Attachment A.

### 3 Methods - Analytical and Preparatory

#### 3.1 Water Sample Preparation for LC/MS/MS

Ten milliliters of sample was transferred into a 50 mL centrifuge tube. Ten milliliters of acetonitrile was added to the sample. After shaking, the sample was sonicated for approximately 2 hours then centrifuged at 3000 rpm for ~10 minutes. A 1 mL portion of the supernatant was transferred to an autosampler vial and fortified with an internal standard solution. The samples were then analyzed using electrospray LC/MS/MS.

#### 3.2 Water Sample Preparation for GC/MS

One hundred milliliters of sample was transferred into a 250 mL polypropylene bottle. Forty milliliter of methyl tert-butyl ether (MTBE) was added to the bottle. The bottle was capped and then shaken for one hour on a reciprocation shaker. The content of the bottle was poured into a 250 mL separatory funnel. The bottle was rinsed with approximately 10 mL fresh MTBE and the rinsate was added to the separatory funnel. The aqueous phase in the funnel was discarded. The organic phase was collected in a 300 mL flask and dried with sodium sulfate. The dried organic phase was then quantitatively transferred into a 50 mL polypropylene centrifuge tube and concentrated to 1 mL using a nitrogen evaporator. The extracted was transferred into a 2 mL GC vial and 10  $\mu$ L of internal standard was added. The sample was then analyzed by GC/MS.

#### 3.3 Sample Analysis by LC/MS/MS

In High Pressure Liquid Chromatography (HPLC), an aliquot of extract is injected and passed through a liquid-phase chromatographic column. Based on the affinity of the analyte for the stationary phase in the column relative to the liquid mobile phase, the analyte is retained for a characteristic amount of time. Following HPLC separation, mass spectrometry provides a rapid and accurate means for analyzing a wide range of organic compounds. Molecules are ionized, fragmented, and detected. The ions characteristic of the compounds are observed and quantitated against calibration standards.

An HP1100 system interfaced to an Applied Biosystems API 4000 and 5000 LC/MS/MS was used to analyze the sample extracts for quantitation. A gradient elution through a Phenomenex Luna  $\mu$  C8(2) Mercury, 20 x 4.0 mm column was used for separation.

The following gradient was performed for C4-C14 acids, PFBS, PFHS, PFHpS, PFOS, PFDS, FOSA, and  $^{13}$ C PFOA (m+4):

|                   |                               |    |  |
|-------------------|-------------------------------|----|--|
| Mobile Phase (A): | 2mM Ammonium Acetate in Water |    |  |
| Mobile Phase (B): | Methanol                      |    |  |
| Time              | %A                            | %B |  |
| 0.0               | 90                            | 10 |  |
| 0.5               | 90                            | 10 |  |
| 2.0               | 10                            | 90 |  |
| 5.0               | 10                            | 90 |  |

|      |    |     |
|------|----|-----|
| 5.1  | 0  | 100 |
| 6.0  | 0  | 100 |
| 6.1  | 90 | 10  |
| 10.0 | 90 | 10  |

The following gradient was performed for MeFOSAA, EtFOSAA, MeFOSE and EtFOSE:

Mobile Phase (A): 2mM Ammonium Acetate in Water  
Mobile Phase (B): Methanol

| <u>Time</u> | <u>%A</u> | <u>%B</u> |
|-------------|-----------|-----------|
| 0.0         | 75        | 25        |
| 0.5         | 75        | 25        |
| 2.0         | 10        | 90        |
| 5.0         | 10        | 90        |
| 5.1         | 0         | 100       |
| 6.0         | 0         | 100       |
| 6.1         | 75        | 25        |
| 10.0        | 75        | 25        |

The following parameters were used for operation of the mass spectrometer:

| <b>Parameter</b>      | <b>Setting</b>   |
|-----------------------|--|
| Ionization Mode       | Electrospray   |
| Polarity              | Negative   |
| Transitions Monitored | 213→169 (C4 Acid)<br>263→219 (C5 Acid)<br>313→269 (C6 Acid)<br>363→319 (C7 Acid)<br>413→369 (C8 Acid)<br>463→419 (C9 Acid)<br>513→469 (C10 Acid)<br>563→519 (C11 Acid)<br>613→569 (C12 Acid)<br>663→619 (C13 Acid)<br>713→669 (C14 Acid)<br>299→80 (PFBS)<br>399→80 (PFHS)<br>449→99 (PFHpS)<br>499→80 (PFOS)<br>599→99 (PFDS)<br>498→78 (FOSA)<br>217→172 (Internal Std. <sup>13</sup> C PFBA (m+4))<br>415→370 (Internal Std. <sup>13</sup> C PFOA (m+2))<br>515→470 (Internal Std. <sup>13</sup> C PFDA (m+2))<br>503→80 (Internal Std. <sup>13</sup> C PFOS (m+4))<br>417→372 (Surrogate <sup>13</sup> C PFOA (m+4))<br>570→419 (MeFOSAA)<br>584→419 (EtFOSAA)<br>616→59 (MeFOSE)<br>630→59 (EtFOSE) |
| Gas Temperature       | 400°C  |

### 3.4 Sample Analysis by GC/MS

The extracts were injected into a gas chromatograph (GC) equipped with a narrow bore capillary column and mass selective detector. The GC was temperature programmed to separate the analytes, and the analytes eluted from the column were introduced to the mass selective detector and identified by comparing retention times and abundances of characteristic masses to that of known standards. Sample concentration was calculated by comparing the response of the characteristic mass relative to that of the calibration curve.

The GC/MS system was operated using the following conditions:

|                             |   |
|-----------------------------|---|
| Instrument                  | Hewlett-Packard model 6890 Series Gas Chromatograph/model 5973 Mass Selective Detector  |
| Column                      | HP-1, 30 m x 0.25 mm ID, 1.00 $\mu$ m df  |
| Oven Temperature            | Hold at 60°C for 4 min., ramp at 20°C/min. to 140°C, ramp at 40°C to 240°C, hold for 5 minutes  |
| Injector Temperature        | 200°C   |
| Transfer Line Temperature   | 280°C   |
| Carrier Gas                 | Helium  |
| Column Flow                 | 1.0 mL/min, Constant  |
| Injection Mode              | Pulsed Splitless, 30psi for 1.5 min.  |
| Injection Liner             | 4 mm ID Single Gooseneck packed with glass wool   |
| Injection Purge Delay       | 1.5 min.  |
| Purge Flow to Split Vent    | 50 mL/min.  |
| Injection Volume            | 2 $\mu$ L   |
| Electron Multiplier Voltage | From ATUNE + 306V   |
| MS Acquisition Mode         | SIM   |
| Ions Monitored              | MFOET (Internal Standard): m/z 448, m/z 466<br>8:1 FTOH (Surrogate): m/z 363, m/z 431<br>6:2 FTOH: m/z 344, m/z 363<br>7:2s FTOH: m/z 319, m/z 355<br>8:2 FTOH: m/z 405, m/z 463<br>10:2 FTOH: m/z 505, m/z 544 |
| Dwell Time                  | 50ms for each ion   |
| MS Temperature              | Quad: 150°C, Source: 230°C  |

## 4 Analysis by LCMSMS

### 4.1 Calibration

A 9-point calibration curve was analyzed at the beginning of the analytical sequence for C4-C14 acids, PFBS, PFHS, PFHpS, PFOS, PFDS, FOSA, and <sup>13</sup>C PFOA (m+4). A continuing calibration verification (CCV) standard (0.250 ng/mL) was used to verify the accuracy of the calibration curve for the duration of the analytical run. At the minimum every tenth sample was a CCV, not including solvent blanks. The calibration curve and the last passing CCV (70-130%) will then bracket acceptable samples. The calibration points were prepared at 0.0125, 0.025, 0.050, 0.100, 0.250, 0.500, 1.0, 2.5 and 5.0 ng/mL (ppb) for LC/MS/MS analysis. The ratio of the analyte concentration to the IS concentration versus the ratio of the analyte instrument response (area) to the IS response (area) was plotted for each point. Using linear

regression with 1/x weighting, the slope, y-intercept and coefficient of determination ( $r^2$ ) were determined. A calibration curve is acceptable if  $r^2 \geq 0.985$ .

A 9-point calibration curve was analyzed at the beginning of the analytical sequence for MeFOSAA, EtFOSSA, MeFOSE and EtFOSE. A continuing calibration verification (CCV) standard (0.250 ng/mL) was used to verify the accuracy of the calibration curve for the duration of the analytical run. At the minimum every tenth sample was a CCV, not including solvent blanks. The calibration curve and the last passing CCV (70-130%) will then bracket acceptable samples. The calibration points were prepared at 0.0125, 0.025, 0.050, 0.100, 0.250, 0.500, 1.0, 2.5 and 5.0 ng/mL (ppb) for LC/MS/MS analysis. The instrument response versus the concentration was plotted for each point. Using quadratic regression with 1/x weighting, the X variable 1 (a), X variable 2 (b), intercept (c) and coefficient of determination ( $r^2$ ) were determined. A calibration curve is acceptable if  $r^2 \geq 0.985$ .

For the results reported here, calibration criteria were met. The calibration curve is included in the raw data in Attachment C.

## 4.2 Surrogates and Internal Standards

$^{13}\text{C}$  labeled-perfluorooctanoic acid ( $^{13}\text{C}$  PFOA (m+4)) is used as a surrogate for the water samples.

$^{13}\text{C}$  PFOA (m+4) recoveries can be found in Attachment B.

$^{13}\text{C}$  PFBA (m+4) is used as the internal standard for the water samples for C4 – C6 Acids.

$^{13}\text{C}$  PFOA (m+2) is used as the internal standard for the water samples for C7 – C9 Acids.

$^{13}\text{C}$  PFDA (m+2) is used as the internal standard for the water samples for C10 – C14 Acids.

$^{13}\text{C}$  PFOS (m+4) is used as the internal standard for the water samples for PFBS, PFHS, PFOS and FOSA

## 4.3 Laboratory Control Spikes

Laboratory control spikes in the analytical set were prepared during each extraction set by adding a known concentration of the analyte to laboratory reagents. Laboratory control spikes are used to assess method accuracy. The laboratory control spikes must show recoveries between 70-130% or the data is rejected. For the results reported here, the laboratory control spikes were within the acceptable range. Laboratory control spike recoveries are given in Attachment B.

## 4.4 Matrix Spikes

Two matrix spikes were prepared for the water samples, one for C4-C14 acids, PFBS, PFHS, PFHpS, PFOS, PFDS, FOSA, and  $^{13}\text{C}$  PFOA (m+4) analysis, and one for MeFOSAA, EtFOSSA, MeFOSE and EtFOSE analysis. They were prepared by adding a known concentration of the target analyte to a separate sample. Matrix spikes are used to assess method accuracy in the matrix. The matrix spike should show recoveries between 70-130%. For the results reported here, the matrix spike was within the acceptable range with the exception of:

L19874-1 (25-1) Spk C at 0.1 ng/mL for EtFOSAA, MeFOSE, EtFOSE and 1.0 ng/mL for C9-C14 acids, PFBS, PFHS, PFOS, FOSA, PFHpS, and PFDS were outside the acceptable

recoveries of 70-130%. Samples were re-extracted. Matrix effect was determined with the following still outside the acceptable recoveries of 70-130%; L19874-1 (25-1) Spk C at 1.0 ng/mL for EtFOSAA, C9-C14 acids, PFBS, PFHS, PFOS, FOSA, PFHpS, and PFDS. L19874-28 (34-1) Spk D was run to determine matrix effect for C9-C14 acids, PFBS, PFHS, PFOS, FOSA, PFHpS, and PFDS. L19874-28 (34-1) passed with acceptable recoveries.

Matrix spike recoveries are given in Attachment B.

#### **4.5 Duplicate**

Laboratory duplicates were not performed as part of this study.

### **5 Analysis by GC/MS**

#### **5.1 System Suitability and Calibration**

Three system suitability standards were analyzed at the beginning of the analytical sequence. The %RSD of the peak area of each analyte should be  $\leq 20$ .

A 6-point calibration was analyzed. The calibration standard analyses were interspersed throughout the analytical sequence. The calibration points were prepared at 0.1, 0.2, 0.5, 1.0, 2.0 and 5.0  $\mu\text{g/mL}$ , which are equivalent to 1, 2, 5, 10, 20 and 50  $\mu\text{g/L}$  (ppb) in samples. A calibration curve is acceptable if  $r^2 \geq 0.985$ .

For the results reported here, system suitability and calibration criteria were met. The system suitability and calibration curve are included in the raw data in Attachment D.

#### **5.2 Surrogate and Internal Standard**

1H,1H-Perfluoro-1-nonanol (8:1 FTOH) was used as surrogate standard. The recoveries of 8:1 FTOH can be found in Attachment B.

2-Perfluorooctyl-[1,1-2H2]-[1,2-13C2]-ethanol (MFOET) was used as internal standard.

#### **5.3 Laboratory Control Spikes**

Laboratory control spikes in the analytical set were prepared during each extraction set by adding a known concentration of the analyte to laboratory reagents. Laboratory control spikes are used to assess method accuracy. The laboratory control spikes must show recoveries between 50-120%. For the results reported here, the laboratory control spikes were within the acceptable range. Laboratory control spike recoveries are given in Attachment B.

#### **5.4 Duplicate**

A laboratory duplicate sample was performed for sample 35-3.

## 5.5 Matrix Spike

Matrix Spike was prepared for sample 35-3 by adding a known concentration of the target analyte to a separate aliquot of sample. Matrix spike is used to assess method accuracy in the matrix. The matrix spike should show recoveries between 50-120%. For the results reported here, the matrix spike was within the acceptable range.

Matrix spike recoveries are given in Attachment B.

## 6 Data Summary

### 6.1 GC/MS Sample Results

Four samples (27-3, 28-3, 32-3, and 35-3) were re-extracted due to surrogate recoveries lower than 50%. Upon re-extraction the surrogate recoveries of the four samples were still below 50%. The low recovery may be due to the matrix effect of the sample. For each of the four samples that were re-extracted, the extraction with the higher recovery was reported.

The laboratory control sample extracted on 1/12/10 was initially analyzed with dataset 011210. The data was rejected due to a mis-injection. This LCS was re-analyzed with dataset 011310 with acceptable results.

The results are reported in parts per billion (ng/mL) on an as-received basis.

Please see Attachment B for a detailed listing of the analytical results.

## 7 Data/Sample Retention

Samples are disposed of one month after the report is issued unless otherwise specified. All electronic data is archived on retrievable media and hard copy reports are stored in data folders maintained by MPI Research. Hardcopy data is stored for a minimum of five years. The client will be notified 30 days prior to the disposal of hardcopy data.

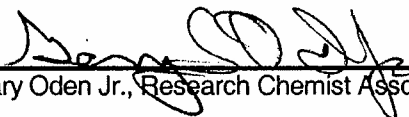
## 8 Attachments

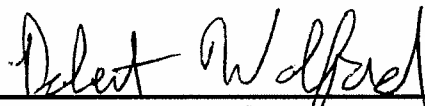
- 7.1 Attachment A: Chain of Custody
- 7.2 Attachment B: Analytical Results
- 7.3 Attachment C: Raw Analytical Data (LC/MS/MS)

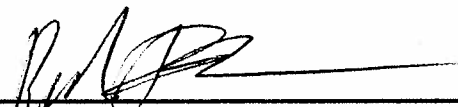


7.4 Attachment D: Raw Analytical Data (GC/MS)

**9 Signatures**

  
\_\_\_\_\_  
Gary Oden Jr., Research Chemist Associate II  
2/2/10  
Date

  
\_\_\_\_\_  
Robert Wolford, Research Chemist Associate I  
2/2/10  
Date

  
\_\_\_\_\_  
Xiang Zhu, Manager, Analytical  
2/2/10  
Date

B

## Fluorochemical Residues in Water Samples By LC/MS/MS

Sample ID: L0019874-1; 25-1

| Analyte  | Amount Found<br>(ng/mL) | LOD (ng/mL) | Analysis Date |
|--|-------------------------|-------------|---------------|
| C4 Acid- Perfluorobutyric Acid                             | ND                      | 0.0125      | 1/25/2010     |
| C5 Acid- Perfluoropentanoic Acid                           | 0.181                   | 0.0125      | 1/25/2010     |
| C6 Acid- Perfluorohexanoic Acid                            | ND                      | 0.0125      | 1/25/2010     |
| C7 Acid- Perfluoroheptanoic Acid                           | ND                      | 0.0125      | 1/25/2010     |
| C8 Acid- Perfluorooctanoic Acid                            | ND                      | 0.0125      | 1/25/2010     |
| C9 Acid- Perfluorononanoic Acid                            | ND                      | 0.0125      | 1/25/2010     |
| C10 Acid- Perfluorodecanoic Acid                           | ND                      | 0.0125      | 1/25/2010     |
| C11 Acid- Perfluoroundecanoic Acid                         | ND                      | 0.0125      | 1/25/2010     |
| C12 Acid- Perfluorododecanoic Acid                         | ND                      | 0.0125      | 1/25/2010     |
| C13 Acid- Perfluorotridecanoic Acid                        | ND                      | 0.0125      | 1/25/2010     |
| C14 Acid- Perfluorotetradecanoic Acid                      | ND                      | 0.0125      | 1/25/2010     |
| PFBS- Perfluorobutanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| PFHS- Perfluorohexanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| PFOS- Perfluorooctanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| FOSA- Perfluorooctane sulfonamide                          | ND                      | 0.0125      | 1/25/2010     |
| PFHpS- Perfluoroheptanesulfonate                           | ND                      | 0.0125      | 1/26/2010     |
| PFDS- Perfluorodecanesulfonate                             | ND                      | 0.0125      | 1/26/2010     |
| MeFOSAA- 2(N-Methylperfluorooctanesulfonamido) acetic acid | ND                      | 0.0125      | 1/25/2010     |
| EtFOSAA- 2(N-Ethylperfluorooctanesulfonamido) acetic acid  | ND                      | 0.0125      | 1/25/2010     |
| MeFOSE- N-Methylperfluorooctane sulfonamidoethanol         | 2.60                    | 0.0125      | 1/29/2010     |
| EtFOSE- N-Ethylperfluorooctane sulfonamidoethanol          | ND                      | 0.0125      | 1/25/2010     |

ND = Not detected = Response is below the LOD of 0.0125 ng/mL (ppb).

NQ = Not quantifiable = Response is between the LOD and the LOQ of 0.0250 ng/mL (ppb).

---

**Fluorotelomer Analysis by GC/MS****Sample ID: L0019874-3; 25-3****Date of Extraction: 1/12/2010****Date Analyzed: 1/12/2010**

| <b>Analyte</b> | <b>Amount Found (ng/mL)</b> | <b>LOQ (ng/mL)</b> |
|----------------|-----------------------------|--------------------|
| 7-2s FTOH      | ND                          | 1.00               |
| 6-2 FTOH       | ND                          | 1.00               |
| 8-2 FTOH       | ND                          | 1.00               |
| 10-2 FTOH      | ND                          | 1.00               |



### Fluorochemical Residues in Water Samples By LC/MS/MS

Sample ID: L0019874-4; 26-1

| Analyte  | Amount Found<br>(ng/mL) | LOD (ng/mL) | Analysis Date |
|--|-------------------------|-------------|---------------|
| C4 Acid- Perfluorobutyric Acid                             | ND                      | 0.0125      | 1/25/2010     |
| C5 Acid- Perfluoropentanoic Acid                           | 0.0483                  | 0.0125      | 1/25/2010     |
| C6 Acid- Perfluorohexanoic Acid                            | 0.0760                  | 0.0125      | 1/25/2010     |
| C7 Acid- Perfluoroheptanoic Acid                           | 0.0370                  | 0.0125      | 1/25/2010     |
| C8 Acid- Perfluorooctanoic Acid                            | 0.142                   | 0.0125      | 1/25/2010     |
| C9 Acid- Perfluorononanoic Acid                            | NQ                      | 0.0125      | 1/25/2010     |
| C10 Acid- Perfluorodecanoic Acid                           | 0.0323                  | 0.0125      | 1/25/2010     |
| C11 Acid- Perfluoroundecanoic Acid                         | ND                      | 0.0125      | 1/25/2010     |
| C12 Acid- Perfluorododecanoic Acid                         | ND                      | 0.0125      | 1/25/2010     |
| C13 Acid- Perfluorotridecanoic Acid                        | ND                      | 0.0125      | 1/25/2010     |
| C14 Acid- Perfluorotetradecanoic Acid                      | ND                      | 0.0125      | 1/25/2010     |
| PFBS- Perfluorobutanesulfonate                             | 0.645                   | 0.0125      | 1/25/2010     |
| PFHS- Perfluorohexanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| PFOS- Perfluorooctanesulfonate                             | 0.538                   | 0.0125      | 1/25/2010     |
| FOSA- Perfluorooctane sulfonamide                          | 0.0430                  | 0.0125      | 1/25/2010     |
| PFHpS- Perfluoroheptanesulfonate                           | ND                      | 0.0125      | 1/25/2010     |
| PFDS- Perfluorodecanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| MeFOSAA- 2(N-Methylperfluorooctanesulfonamido) acetic acid | 0.0686                  | 0.0125      | 1/25/2010     |
| EtFOSAA- 2(N-Ethylperfluorooctanesulfonamido) acetic acid  | 0.0460                  | 0.0125      | 1/25/2010     |
| MeFOSE- N-Methylperfluorooctane sulfonamidoethanol         | 3.65                    | 0.0125      | 1/25/2010     |
| EtFOSE- N-Ethylperfluorooctane sulfonamidoethanol          | ND                      | 0.0125      | 1/25/2010     |

ND = Not detected = Response is below the LOD of 0.0125 ng/mL (ppb).

NQ = Not quantifiable = Response is between the LOD and the LOQ of 0.0250 ng/mL (ppb).



3058 Research Drive  
State College, Pennsylvania 16801 USA  
Telephone: 814.272.1039  
Fax: 814.272.1019

## Analytical Report

---

### Fluorotelomer Analysis by GC/MS

**Sample ID: L0019847-6; 26-3**

**Date of Extraction: 1/12/2010**

**Date Analyzed: 1/12/2010**

| Analyte   | Amount Found (ng/mL) | LOQ (ng/mL) |
|-----------|----------------------|-------------|
| 7-2s FTOH | ND                   | 1.00        |
| 6-2 FTOH  | ND                   | 1.00        |
| 8-2 FTOH  | ND                   | 1.00        |
| 10-2 FTOH | ND                   | 1.00        |

## Fluorochemical Residues in Water Samples By LC/MS/MS

Sample ID: L0019874-7; 27-1

| Analyte  | Amount Found<br>(ng/mL) | LOD (ng/mL) | Analysis Date |
|--|-------------------------|-------------|---------------|
| C4 Acid- Perfluorobutyric Acid                             | ND                      | 0.0125      | 1/25/2010     |
| C5 Acid- Perfluoropentanoic Acid                           | 0.0528                  | 0.0125      | 1/25/2010     |
| C6 Acid- Perfluorohexanoic Acid                            | 0.0598                  | 0.0125      | 1/25/2010     |
| C7 Acid- Perfluoroheptanoic Acid                           | 0.0532                  | 0.0125      | 1/25/2010     |
| C8 Acid- Perfluorooctanoic Acid                            | 0.134                   | 0.0125      | 1/25/2010     |
| C9 Acid- Perfluorononanoic Acid                            | 0.0531                  | 0.0125      | 1/25/2010     |
| C10 Acid- Perfluorodecanoic Acid                           | NQ                      | 0.0125      | 1/25/2010     |
| C11 Acid- Perfluoroundecanoic Acid                         | ND                      | 0.0125      | 1/25/2010     |
| C12 Acid- Perfluorododecanoic Acid                         | ND                      | 0.0125      | 1/25/2010     |
| C13 Acid- Perfluorotridecanoic Acid                        | ND                      | 0.0125      | 1/25/2010     |
| C14 Acid- Perfluorotetradecanoic Acid                      | ND                      | 0.0125      | 1/25/2010     |
| PFBS- Perfluorobutanesulfonate                             | 0.0486                  | 0.0125      | 1/25/2010     |
| PFHS- Perfluorohexanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| PFOS- Perfluorooctanesulfonate                             | 1.18                    | 0.0125      | 1/25/2010     |
| FOSA- Perfluorooctane sulfonamide                          | ND                      | 0.0125      | 1/25/2010     |
| PFHpS- Perfluoroheptanesulfonate                           | ND                      | 0.0125      | 1/25/2010     |
| PFDS- Perfluorodecanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| MeFOSAA- 2(N-Methylperfluorooctanesulfonamido) acetic acid | 0.123                   | 0.0125      | 1/25/2010     |
| EtFOSAA- 2(N-Ethylperfluorooctanesulfonamido) acetic acid  | ND                      | 0.0125      | 1/25/2010     |
| MeFOSE- N-Methylperfluorooctane sulfonamidoethanol         | 0.137                   | 0.0125      | 1/25/2010     |
| EtFOSE- N-Ethylperfluorooctane sulfonamidoethanol          | ND                      | 0.0125      | 1/25/2010     |

ND = Not detected = Response is below the LOD of 0.0125 ng/mL (ppb).

NQ = Not quantifiable = Response is between the LOD and the LOQ of 0.0250 ng/mL (ppb).



3058 Research Drive  
State College, Pennsylvania 16801 USA  
Telephone: 814.272.1039  
Fax: 814.272.1019

## Analytical Report

---

### Fluorotelomer Analysis by GC/MS

**Sample ID: L0019874-9; 27-3**

**Date of Extraction: 1/19/2010**

**Date Analyzed: 1/19/2010**

| Analyte   | Amount Found (ng/mL) | LOQ (ng/mL) |
|-----------|----------------------|-------------|
| 7-2s FTOH | ND                   | 1.00        |
| 6-2 FTOH  | ND                   | 1.00        |
| 8-2 FTOH  | ND                   | 1.00        |
| 10-2 FTOH | ND                   | 1.00        |



## Fluorochemical Residues in Water Samples By LC/MS/MS

Sample ID: L0019874-10; 28-1

| Analyte  | Amount Found<br>(ng/mL) | LOD (ng/mL) | Analysis Date |
|--|-------------------------|-------------|---------------|
| C4 Acid- Perfluorobutyric Acid                             | ND                      | 0.0125      | 1/25/2010     |
| C5 Acid- Perfluoropentanoic Acid                           | 0.0438                  | 0.0125      | 1/25/2010     |
| C6 Acid- Perfluorohexanoic Acid                            | 0.0604                  | 0.0125      | 1/25/2010     |
| C7 Acid- Perfluoroheptanoic Acid                           | 0.0582#                 | 0.0125      | 1/25/2010     |
| C8 Acid- Perfluorooctanoic Acid                            | 0.124#                  | 0.0125      | 1/25/2010     |
| C9 Acid- Perfluorononanoic Acid                            | 0.163#                  | 0.0125      | 1/25/2010     |
| C10 Acid- Perfluorodecanoic Acid                           | 0.0365                  | 0.0125      | 1/25/2010     |
| C11 Acid- Perfluoroundecanoic Acid                         | NQ                      | 0.0125      | 1/25/2010     |
| C12 Acid- Perfluorododecanoic Acid                         | ND                      | 0.0125      | 1/25/2010     |
| C13 Acid- Perfluorotridecanoic Acid                        | ND                      | 0.0125      | 1/25/2010     |
| C14 Acid- Perfluorotetradecanoic Acid                      | ND                      | 0.0125      | 1/25/2010     |
| PFBS- Perfluorobutanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| PFHS- Perfluorohexanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| PFOS- Perfluorooctanesulfonate                             | 24.2                    | 0.0125      | 1/25/2010     |
| FOSA- Perfluorooctane sulfonamide                          | ND                      | 0.0125      | 1/25/2010     |
| PFHpS- Perfluoroheptanesulfonate                           | ND                      | 0.0125      | 1/25/2010     |
| PFDS- Perfluorodecanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| MeFOSAA- 2(N-Methylperfluorooctanesulfonamido) acetic acid | 0.0301                  | 0.0125      | 1/25/2010     |
| EtFOSAA- 2(N-Ethylperfluorooctanesulfonamido) acetic acid  | ND                      | 0.0125      | 1/25/2010     |
| MeFOSE- N-Methylperfluorooctane sulfonamidoethanol         | ND                      | 0.0125      | 1/25/2010     |
| EtFOSE- N-Ethylperfluorooctane sulfonamidoethanol          | ND                      | 0.0125      | 1/25/2010     |

ND = Not detected = Response is below the LOD of 0.0125 ng/mL (ppb).

NQ = Not quantifiable = Response is between the LOD and the LOQ of 0.0250 ng/mL (ppb).



3058 Research Drive  
State College, Pennsylvania 16801 USA  
Telephone: 814.272.1039  
Fax: 814.272.1019

## Analytical Report

---

### Fluorotelomer Analysis by GC/MS

**Sample ID: L0019874-12; 28-3**

**Date of Extraction: 1/19/2010**

**Date Analyzed: 1/19/2010**

| Analyte   | Amount Found (ng/mL) | LOQ (ng/mL) |
|-----------|----------------------|-------------|
| 7-2s FTOH | ND                   | 1.00        |
| 6-2 FTOH  | ND                   | 1.00        |
| 8-2 FTOH  | 1.40                 | 1.00        |
| 10-2 FTOH | ND                   | 1.00        |

## Fluorochemical Residues in Water Samples By LC/MS/MS

Sample ID: L0019874-13; 29-1

| Analyte  | Amount Found<br>(ng/mL) | LOD (ng/mL) | Analysis Date |
|--|-------------------------|-------------|---------------|
| C4 Acid- Perfluorobutyric Acid                             | ND                      | 0.0125      | 1/25/2010     |
| C5 Acid- Perfluoropentanoic Acid                           | ND                      | 0.0125      | 1/25/2010     |
| C6 Acid- Perfluorohexanoic Acid                            | NQ                      | 0.0125      | 1/25/2010     |
| C7 Acid- Perfluoroheptanoic Acid                           | ND                      | 0.0125      | 1/25/2010     |
| C8 Acid- Perfluorooctanoic Acid                            | 0.0445                  | 0.0125      | 1/25/2010     |
| C9 Acid- Perfluorononanoic Acid                            | ND                      | 0.0125      | 1/25/2010     |
| C10 Acid- Perfluorodecanoic Acid                           | ND                      | 0.0125      | 1/25/2010     |
| C11 Acid- Perfluoroundecanoic Acid                         | ND                      | 0.0125      | 1/25/2010     |
| C12 Acid- Perfluorododecanoic Acid                         | ND                      | 0.0125      | 1/25/2010     |
| C13 Acid- Perfluorotridecanoic Acid                        | ND                      | 0.0125      | 1/25/2010     |
| C14 Acid- Perfluorotetradecanoic Acid                      | ND                      | 0.0125      | 1/25/2010     |
| PFBS- Perfluorobutanesulfonate                             | 0.0569                  | 0.0125      | 1/25/2010     |
| PFHS- Perfluorohexanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| PFOS- Perfluorooctanesulfonate                             | 0.234                   | 0.0125      | 1/25/2010     |
| FOSA- Perfluorooctane sulfonamide                          | ND                      | 0.0125      | 1/25/2010     |
| PFHpS- Perfluoroheptanesulfonate                           | ND                      | 0.0125      | 1/25/2010     |
| PFDS- Perfluorodecanesulfonate                             | 0.0563                  | 0.0125      | 1/25/2010     |
| MeFOSAA- 2(N-Methylperfluorooctanesulfonamido) acetic acid | ND                      | 0.0125      | 1/25/2010     |
| EtFOSAA- 2(N-Ethylperfluorooctanesulfonamido) acetic acid  | ND                      | 0.0125      | 1/25/2010     |
| MeFOSE- N-Methylperfluorooctane sulfonamidoethanol         | 0.124                   | 0.0125      | 1/25/2010     |
| EtFOSE- N-Ethylperfluorooctane sulfonamidoethanol          | 0.0430                  | 0.0125      | 1/25/2010     |

ND = Not detected = Response is below the LOD of 0.0125 ng/mL (ppb).

NQ = Not quantifiable = Response is between the LOD and the LOQ of 0.0250 ng/mL (ppb).

---

**Fluorotelomer Analysis by GC/MS****Sample ID: L0019874-15; 29-3****Date of Extraction: 1/12/2010****Date Analyzed: 1/12/2010**

| Analyte   | Amount Found (ng/mL) | LOQ (ng/mL) |
|-----------|----------------------|-------------|
| 7-2s FTOH | ND                   | 1.00        |
| 6-2 FTOH  | 1.40                 | 1.00        |
| 8-2 FTOH  | ND                   | 1.00        |
| 10-2 FTOH | ND                   | 1.00        |

## Fluorochemical Residues in Water Samples By LC/MS/MS

Sample ID: L0019874-16; 30-1

| Analyte  | Amount Found (ng/mL) | LOD (ng/mL) | Analysis Date |
|--|----------------------|-------------|---------------|
| C4 Acid- Perfluorobutyric Acid                             | ND                   | 0.0125      | 1/25/2010     |
| C5 Acid- Perfluoropentanoic Acid                           | 0.160                | 0.0125      | 1/25/2010     |
| C6 Acid- Perfluorohexanoic Acid                            | 0.130                | 0.0125      | 1/25/2010     |
| C7 Acid- Perfluoroheptanoic Acid                           | 0.0987               | 0.0125      | 1/25/2010     |
| C8 Acid- Perfluorooctanoic Acid                            | 0.294                | 0.0125      | 1/25/2010     |
| C9 Acid- Perfluorononanoic Acid                            | 0.0454               | 0.0125      | 1/25/2010     |
| C10 Acid- Perfluorodecanoic Acid                           | 0.0679               | 0.0125      | 1/25/2010     |
| C11 Acid- Perfluoroundecanoic Acid                         | NQ                   | 0.0125      | 1/25/2010     |
| C12 Acid- Perfluorododecanoic Acid                         | ND                   | 0.0125      | 1/25/2010     |
| C13 Acid- Perfluorotridecanoic Acid                        | ND                   | 0.0125      | 1/25/2010     |
| C14 Acid- Perfluorotetradecanoic Acid                      | ND                   | 0.0125      | 1/25/2010     |
| PFBS- Perfluorobutanesulfonate                             | NQ                   | 0.0125      | 1/25/2010     |
| PFHS- Perfluorohexanesulfonate                             | ND                   | 0.0125      | 1/25/2010     |
| PFOS- Perfluorooctanesulfonate                             | 1.13                 | 0.0125      | 1/25/2010     |
| FOSA- Perfluorooctane sulfonamide                          | 0.0387               | 0.0125      | 1/25/2010     |
| PFHpS- Perfluoroheptanesulfonate                           | ND                   | 0.0125      | 1/25/2010     |
| PFDS- Perfluorodecanesulfonate                             | ND                   | 0.0125      | 1/25/2010     |
| MeFOSAA- 2(N-Methylperfluorooctanesulfonamido) acetic acid | 0.342                | 0.0125      | 1/25/2010     |
| EtFOSAA- 2(N-Ethylperfluorooctanesulfonamido) acetic acid  | 0.0677               | 0.0125      | 1/25/2010     |
| MeFOSE- N-Methylperfluorooctane sulfonamidoethanol         | 0.0401               | 0.0125      | 1/25/2010     |
| EtFOSE- N-Ethylperfluorooctane sulfonamidoethanol          | ND                   | 0.0125      | 1/25/2010     |

ND = Not detected = Response is below the LOD of 0.0125 ng/mL (ppb).

NQ = Not quantifiable = Response is between the LOD and the LOQ of 0.0250 ng/mL (ppb).



3058 Research Drive  
State College, Pennsylvania 16801 USA  
Telephone: 814.272.1039  
Fax: 814.272.1019

## Analytical Report

---

### Fluorotelomer Analysis by GC/MS

**Sample ID: L0019874-18; 30-3**

**Date of Extraction: 1/12/2010**

**Date Analyzed: 1/12/2010**

| Analyte   | Amount Found (ng/mL) | LOQ (ng/mL) |
|-----------|----------------------|-------------|
| 7-2s FTOH | ND                   | 1.00        |
| 6-2 FTOH  | 2.00                 | 1.00        |
| 8-2 FTOH  | ND                   | 1.00        |
| 10-2 FTOH | ND                   | 1.00        |

## Fluorochemical Residues in Water Samples By LC/MS/MS

Sample ID: L0019874-19; 31-1

| Analyte  | Amount Found<br>(ng/mL) | LOD (ng/mL) | Analysis Date |
|--|-------------------------|-------------|---------------|
| C4 Acid- Perfluorobutyric Acid                             | ND                      | 0.0125      | 1/25/2010     |
| C5 Acid- Perfluoropentanoic Acid                           | ND                      | 0.0125      | 1/25/2010     |
| C6 Acid- Perfluorohexanoic Acid                            | ND                      | 0.0125      | 1/25/2010     |
| C7 Acid- Perfluoroheptanoic Acid                           | ND                      | 0.0125      | 1/25/2010     |
| C8 Acid- Perfluorooctanoic Acid                            | ND                      | 0.0125      | 1/25/2010     |
| C9 Acid- Perfluorononanoic Acid                            | 0.0391                  | 0.0125      | 1/25/2010     |
| C10 Acid- Perfluorodecanoic Acid                           | ND                      | 0.0125      | 1/25/2010     |
| C11 Acid- Perfluoroundecanoic Acid                         | ND                      | 0.0125      | 1/25/2010     |
| C12 Acid- Perfluorododecanoic Acid                         | ND                      | 0.0125      | 1/25/2010     |
| C13 Acid- Perfluorotridecanoic Acid                        | ND                      | 0.0125      | 1/25/2010     |
| C14 Acid- Perfluorotetradecanoic Acid                      | ND                      | 0.0125      | 1/25/2010     |
| PFBS- Perfluorobutanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| PFHS- Perfluorohexanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| PFOS- Perfluorooctanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| FOSA- Perfluorooctane sulfonamide                          | ND                      | 0.0125      | 1/25/2010     |
| PFHpS- Perfluoroheptanesulfonate                           | ND                      | 0.0125      | 1/25/2010     |
| PFDS- Perfluorodecanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| MeFOSAA- 2(N-Methylperfluorooctanesulfonamido) acetic acid | ND                      | 0.0125      | 1/25/2010     |
| EtFOSAA- 2(N-Ethylperfluorooctanesulfonamido) acetic acid  | ND                      | 0.0125      | 1/25/2010     |
| MeFOSE- N-Methylperfluorooctane sulfonamidoethanol         | 0.573                   | 0.0125      | 1/25/2010     |
| EtFOSE- N-Ethylperfluorooctane sulfonamidoethanol          | ND                      | 0.0125      | 1/25/2010     |

ND = Not detected = Response is below the LOD of 0.0125 ng/mL (ppb).

NQ = Not quantifiable = Response is between the LOD and the LOQ of 0.0250 ng/mL (ppb).

---

**Fluorotelomer Analysis by GC/MS**

**Sample ID: L0019874-21; 31-3**

**Date of Extraction: 1/12/2010**

**Date Analyzed: 1/12/2010**

| Analyte   | Amount Found (ng/mL) | LOQ (ng/mL) |
|-----------|----------------------|-------------|
| 7-2s FTOH | ND                   | 1.00        |
| 6-2 FTOH  | ND                   | 1.00        |
| 8-2 FTOH  | ND                   | 1.00        |
| 10-2 FTOH | ND                   | 1.00        |





**Fluorochemical Residues in Water Samples By LC/MS/MS**

Sample ID: L0019874-22; 32-1

| Analyte  | Amount Found<br>(ng/mL) | LOD (ng/mL) | Analysis Date |
|--|-------------------------|-------------|---------------|
| C4 Acid- Perfluorobutyric Acid                             | ND                      | 0.0125      | 1/25/2010     |
| C5 Acid- Perfluoropentanoic Acid                           | NQ                      | 0.0125      | 1/25/2010     |
| C6 Acid- Perfluorohexanoic Acid                            | ND                      | 0.0125      | 1/25/2010     |
| C7 Acid- Perfluoroheptanoic Acid                           | ND                      | 0.0125      | 1/25/2010     |
| C8 Acid- Perfluorooctanoic Acid                            | ND                      | 0.0125      | 1/25/2010     |
| C9 Acid- Perfluorononanoic Acid                            | 0.0367                  | 0.0125      | 1/25/2010     |
| C10 Acid- Perfluorodecanoic Acid                           | ND                      | 0.0125      | 1/25/2010     |
| C11 Acid- Perfluoroundecanoic Acid                         | ND                      | 0.0125      | 1/25/2010     |
| C12 Acid- Perfluorododecanoic Acid                         | ND                      | 0.0125      | 1/25/2010     |
| C13 Acid- Perfluorotridecanoic Acid                        | ND                      | 0.0125      | 1/25/2010     |
| C14 Acid- Perfluorotetradecanoic Acid                      | ND                      | 0.0125      | 1/25/2010     |
| PFBS- Perfluorobutanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| PFHS- Perfluorohexanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| PFOS- Perfluorooctanesulfonate                             | 5.99                    | 0.0125      | 1/25/2010     |
| FOSA- Perfluorooctane sulfonamide                          | ND                      | 0.0125      | 1/25/2010     |
| PFHpS- Perfluoroheptanesulfonate                           | ND                      | 0.0125      | 1/25/2010     |
| PFDS- Perfluorodecanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| MeFOSAA- 2(N-Methylperfluorooctanesulfonamido) acetic acid | ND                      | 0.0125      | 1/25/2010     |
| EtFOSAA- 2(N-Ethylperfluorooctanesulfonamido) acetic acid  | ND                      | 0.0125      | 1/25/2010     |
| MeFOSE- N-Methylperfluorooctane sulfonamidoethanol         | ND                      | 0.0125      | 1/25/2010     |
| EtFOSE- N-Ethylperfluorooctane sulfonamidoethanol          | ND                      | 0.0125      | 1/25/2010     |

ND = Not detected = Response is below the LOD of 0.0125 ng/mL (ppb).

NQ = Not quantifiable = Response is between the LOD and the LOQ of 0.0250 ng/mL (ppb).



3058 Research Drive  
State College, Pennsylvania 16801 USA  
Telephone: 814.272.1039  
Fax: 814.272.1019

## Analytical Report

---

### Fluorotelomer Analysis by GC/MS

**Sample ID: L0019874-24; 32-3**

**Date of Extraction: 1/12/2010**

**Date Analyzed: 1/12/2010**

| Analyte   | Amount Found (ng/mL) | LOQ (ng/mL) |
|-----------|----------------------|-------------|
| 7-2s FTOH | ND                   | 1.00        |
| 6-2 FTOH  | ND                   | 1.00        |
| 8-2 FTOH  | ND                   | 1.00        |
| 10-2 FTOH | ND                   | 1.00        |



## Fluorochemical Residues in Water Samples By LC/MS/MS

Sample ID: L0019874-25; 33-1

| Analyte  | Amount Found<br>(ng/mL) | LOD (ng/mL) | Analysis Date |
|--|-------------------------|-------------|---------------|
| C4 Acid- Perfluorobutyric Acid                             | ND                      | 0.0125      | 1/25/2010     |
| C5 Acid- Perfluoropentanoic Acid                           | ND                      | 0.0125      | 1/25/2010     |
| C6 Acid- Perfluorohexanoic Acid                            | NQ                      | 0.0125      | 1/25/2010     |
| C7 Acid- Perfluoroheptanoic Acid                           | ND                      | 0.0125      | 1/25/2010     |
| C8 Acid- Perfluorooctanoic Acid                            | NQ                      | 0.0125      | 1/25/2010     |
| C9 Acid- Perfluorononanoic Acid                            | ND                      | 0.0125      | 1/25/2010     |
| C10 Acid- Perfluorodecanoic Acid                           | ND                      | 0.0125      | 1/25/2010     |
| C11 Acid- Perfluoroundecanoic Acid                         | ND                      | 0.0125      | 1/25/2010     |
| C12 Acid- Perfluorododecanoic Acid                         | ND                      | 0.0125      | 1/25/2010     |
| C13 Acid- Perfluorotridecanoic Acid                        | ND                      | 0.0125      | 1/25/2010     |
| C14 Acid- Perfluorotetradecanoic Acid                      | ND                      | 0.0125      | 1/25/2010     |
| PFBS- Perfluorobutanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| PFHS- Perfluorohexanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| PFOS- Perfluorooctanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| FOSA- Perfluorooctane sulfonamide                          | ND                      | 0.0125      | 1/25/2010     |
| PFHpS- Perfluoroheptanesulfonate                           | ND                      | 0.0125      | 1/25/2010     |
| PFDS- Perfluorodecanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| MeFOSAA- 2(N-Methylperfluorooctanesulfonamido) acetic acid | ND                      | 0.0125      | 1/25/2010     |
| EtFOSAA- 2(N-Ethylperfluorooctanesulfonamido) acetic acid  | ND                      | 0.0125      | 1/25/2010     |
| MeFOSE- N-Methylperfluorooctane sulfonamidoethanol         | ND                      | 0.0125      | 1/25/2010     |
| EtFOSE- N-Ethylperfluorooctane sulfonamidoethanol          | ND                      | 0.0125      | 1/25/2010     |

ND = Not detected = Response is below the LOD of 0.0125 ng/mL (ppb).

NQ = Not quantifiable = Response is between the LOD and the LOQ of 0.0250 ng/mL (ppb).



3058 Research Drive  
State College, Pennsylvania 16801 USA  
Telephone: 814.272.1039  
Fax: 814.272.1019

## Analytical Report

---

### Fluorotelomer Analysis by GC/MS

**Sample ID: L0019874-27; 33-3**

**Date of Extraction: 1/12/2010**

**Date Analyzed: 1/13/2010**

| Analyte   | Amount Found (ng/mL) | LOQ (ng/mL) |
|-----------|----------------------|-------------|
| 7-2s FTOH | ND                   | 1.00        |
| 6-2 FTOH  | ND                   | 1.00        |
| 8-2 FTOH  | ND                   | 1.00        |
| 10-2 FTOH | ND                   | 1.00        |



## Fluorochemical Residues in Water Samples By LC/MS/MS

Sample ID: L0019874-28; 34-1

| Analyte  | Amount Found<br>(ng/mL) | LOD (ng/mL) | Analysis Date |
|--|-------------------------|-------------|---------------|
| C4 Acid- Perfluorobutyric Acid                             | ND                      | 0.0125      | 1/25/2010     |
| C5 Acid- Perfluoropentanoic Acid                           | ND                      | 0.0125      | 1/25/2010     |
| C6 Acid- Perfluorohexanoic Acid                            | 0.0478                  | 0.0125      | 1/25/2010     |
| C7 Acid- Perfluoroheptanoic Acid                           | NQ                      | 0.0125      | 1/25/2010     |
| C8 Acid- Perfluorooctanoic Acid                            | 0.0461                  | 0.0125      | 1/25/2010     |
| C9 Acid- Perfluorononanoic Acid                            | ND                      | 0.0125      | 1/25/2010     |
| C10 Acid- Perfluorodecanoic Acid                           | ND                      | 0.0125      | 1/25/2010     |
| C11 Acid- Perfluoroundecanoic Acid                         | ND                      | 0.0125      | 1/25/2010     |
| C12 Acid- Perfluorododecanoic Acid                         | ND                      | 0.0125      | 1/25/2010     |
| C13 Acid- Perfluorotridecanoic Acid                        | ND                      | 0.0125      | 1/25/2010     |
| C14 Acid- Perfluorotetradecanoic Acid                      | ND                      | 0.0125      | 1/25/2010     |
| PFBS- Perfluorobutanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| PFHS- Perfluorohexanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| PFOS- Perfluorooctanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| FOSA- Perfluorooctane sulfonamide                          | ND                      | 0.0125      | 1/25/2010     |
| PFHpS- Perfluoroheptanesulfonate                           | ND                      | 0.0125      | 1/25/2010     |
| PFDS- Perfluorodecanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| MeFOSAA- 2(N-Methylperfluorooctanesulfonamido) acetic acid | ND                      | 0.0125      | 1/25/2010     |
| EtFOSAA- 2(N-Ethylperfluorooctanesulfonamido) acetic acid  | ND                      | 0.0125      | 1/25/2010     |
| MeFOSE- N-Methylperfluorooctane sulfonamidoethanol         | ND                      | 0.0125      | 1/25/2010     |
| EtFOSE- N-Ethylperfluorooctane sulfonamidoethanol          | ND                      | 0.0125      | 1/25/2010     |

ND = Not detected = Response is below the LOD of 0.0125 ng/mL (ppb).

NQ = Not quantifiable = Response is between the LOD and the LOQ of 0.0250 ng/mL (ppb).



3058 Research Drive  
State College, Pennsylvania 16801 USA  
Telephone: 814.272.1039  
Fax: 814.272.1019

## Analytical Report

---

### Fluorotelomer Analysis by GC/MS

**Sample ID: L0019874-30; 34-3**

**Date of Extraction: 1/12/2010**

**Date Analyzed: 1/13/2010**

| Analyte   | Amount Found (ng/mL) | LOQ (ng/mL) |
|-----------|----------------------|-------------|
| 7-2s FTOH | ND                   | 1.00        |
| 6-2 FTOH  | ND                   | 1.00        |
| 8-2 FTOH  | ND                   | 1.00        |
| 10-2 FTOH | ND                   | 1.00        |



**Fluorochemical Residues in Water Samples By LC/MS/MS**

Sample ID: L0019874-31; 35-1

| Analyte  | Amount Found<br>(ng/mL) | LOD (ng/mL) | Analysis Date |
|--|-------------------------|-------------|---------------|
| C4 Acid- Perfluorobutyric Acid                             | ND                      | 0.0125      | 1/25/2010     |
| C5 Acid- Perfluoropentanoic Acid                           | ND                      | 0.0125      | 1/25/2010     |
| C6 Acid- Perfluorohexanoic Acid                            | 0.0529                  | 0.0125      | 1/25/2010     |
| C7 Acid- Perfluoroheptanoic Acid                           | 0.0533                  | 0.0125      | 1/25/2010     |
| C8 Acid- Perfluorooctanoic Acid                            | 0.0635                  | 0.0125      | 1/25/2010     |
| C9 Acid- Perfluorononanoic Acid                            | 0.0256                  | 0.0125      | 1/25/2010     |
| C10 Acid- Perfluorodecanoic Acid                           | ND                      | 0.0125      | 1/25/2010     |
| C11 Acid- Perfluoroundecanoic Acid                         | ND                      | 0.0125      | 1/25/2010     |
| C12 Acid- Perfluorododecanoic Acid                         | ND                      | 0.0125      | 1/25/2010     |
| C13 Acid- Perfluorotridecanoic Acid                        | ND                      | 0.0125      | 1/25/2010     |
| C14 Acid- Perfluorotetradecanoic Acid                      | ND                      | 0.0125      | 1/25/2010     |
| PFBS- Perfluorobutanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| PFHS- Perfluorohexanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| PFOS- Perfluorooctanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| FOSA- Perfluorooctane sulfonamide                          | ND                      | 0.0125      | 1/25/2010     |
| PFHpS- Perfluoroheptanesulfonate                           | ND                      | 0.0125      | 1/25/2010     |
| PFDS- Perfluorodecanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| MeFOSAA- 2(N-Methylperfluorooctanesulfonamido) acetic acid | ND                      | 0.0125      | 1/25/2010     |
| EtFOSAA- 2(N-Ethylperfluorooctanesulfonamido) acetic acid  | ND                      | 0.0125      | 1/25/2010     |
| MeFOSE- N-Methylperfluorooctane sulfonamidoethanol         | ND                      | 0.0125      | 1/25/2010     |
| EtFOSE- N-Ethylperfluorooctane sulfonamidoethanol          | ND                      | 0.0125      | 1/25/2010     |

ND = Not detected = Response is below the LOD of 0.0125 ng/mL (ppb).

NQ = Not quantifiable = Response is between the LOD and the LOQ of 0.0250 ng/mL (ppb).



3058 Research Drive  
State College, Pennsylvania 16801 USA  
Telephone: 814.272.1039  
Fax: 814.272.1019

## Analytical Report

---

### Fluorotelomer Analysis by GC/MS

**Sample ID: L0019874-33; 35-3**

**Date of Extraction: 1/19/2010**

**Date Analyzed: 1/19/2010**

| Analyte   | Amount Found (ng/mL) | LOQ (ng/mL) |
|-----------|----------------------|-------------|
| 7-2s FTOH | ND                   | 1.00        |
| 6-2 FTOH  | ND                   | 1.00        |
| 8-2 FTOH  | ND                   | 1.00        |
| 10-2 FTOH | ND                   | 1.00        |





**MPI**  
RESEARCH

3058 Research Drive  
State College, Pennsylvania 16801 USA  
Telephone: 814.272.1039  
Fax: 814.272.1019

## Analytical Report

---

### Fluorotelomer Analysis by GC/MS

**Sample ID: L0019874-33 Duplicate; 35-3**

**Date of Extraction: 1/19/2010**

**Date Analyzed: 1/19/2010**

| Analyte   | Amount Found (ng/mL) | LOQ (ng/mL) |
|-----------|----------------------|-------------|
| 7-2s FTOH | ND                   | 1.00        |
| 6-2 FTOH  | ND                   | 1.00        |
| 8-2 FTOH  | ND                   | 1.00        |
| 10-2 FTOH | ND                   | 1.00        |

## Fluorochemical Residues in Water Samples By LC/MS/MS

Sample ID: L0019874-34; 4D1 Duplicate

| Analyte  | Amount Found<br>(ng/mL) | LOD (ng/mL) | Analysis Date |
|--|-------------------------|-------------|---------------|
| C4 Acid- Perfluorobutyric Acid                             | ND                      | 0.0125      | 1/25/2010     |
| C5 Acid- Perfluoropentanoic Acid                           | ND                      | 0.0125      | 1/25/2010     |
| C6 Acid- Perfluorohexanoic Acid                            | ND                      | 0.0125      | 1/25/2010     |
| C7 Acid- Perfluoroheptanoic Acid                           | ND                      | 0.0125      | 1/25/2010     |
| C8 Acid- Perfluorooctanoic Acid                            | ND                      | 0.0125      | 1/25/2010     |
| C9 Acid- Perfluorononanoic Acid                            | 0.0364                  | 0.0125      | 1/25/2010     |
| C10 Acid- Perfluorodecanoic Acid                           | ND                      | 0.0125      | 1/25/2010     |
| C11 Acid- Perfluoroundecanoic Acid                         | ND                      | 0.0125      | 1/25/2010     |
| C12 Acid- Perfluorododecanoic Acid                         | ND                      | 0.0125      | 1/25/2010     |
| C13 Acid- Perfluorotridecanoic Acid                        | ND                      | 0.0125      | 1/25/2010     |
| C14 Acid- Perfluorotetradecanoic Acid                      | ND                      | 0.0125      | 1/25/2010     |
| PFBS- Perfluorobutanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| PFHS- Perfluorohexanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| PFOS- Perfluorooctanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| FOSA- Perfluorooctane sulfonamide                          | ND                      | 0.0125      | 1/25/2010     |
| PFHpS- Perfluoroheptanesulfonate                           | ND                      | 0.0125      | 1/25/2010     |
| PFDS- Perfluorodecanesulfonate                             | ND                      | 0.0125      | 1/25/2010     |
| MeFOSAA- 2(N-Methylperfluorooctanesulfonamido) acetic acid | ND                      | 0.0125      | 1/25/2010     |
| EtFOSAA- 2(N-Ethylperfluorooctanesulfonamido) acetic acid  | ND                      | 0.0125      | 1/25/2010     |
| MeFOSE- N-Methylperfluorooctane sulfonamidoethanol         | 0.460                   | 0.0125      | 1/25/2010     |
| EtFOSE- N-Ethylperfluorooctane sulfonamidoethanol          | ND                      | 0.0125      | 1/25/2010     |

ND = Not detected = Response is below the LOD of 0.0125 ng/mL (ppb).

NQ = Not quantifiable = Response is between the LOD and the LOQ of 0.0250 ng/mL (ppb).



**MPI**  
RESEARCH

3058 Research Drive  
State College, Pennsylvania 16801 USA  
Telephone: 814.272.1039  
Fax: 814.272.1019

## Analytical Report

---

### Fluorotelomer Analysis by GC/MS

**Sample ID: L0019874-36; 4D3 Duplicate**

**Date of Extraction: 1/12/2010**

**Date Analyzed: 1/12/2010**

| Analyte   | Amount Found (ng/mL) | LOQ (ng/mL) |
|-----------|----------------------|-------------|
| 7-2s FTOH | ND                   | 1.00        |
| 6-2 FTOH  | ND                   | 1.00        |
| 8-2 FTOH  | ND                   | 1.00        |
| 10-2 FTOH | ND                   | 1.00        |

## Recovery Summary of Fluorochemical Residues in Water Samples

| Sample Description   | Amount Spiked (ng/mL) | C4 Acid                     |                          |              | C5 Acid                     |                          |              | C6 Acid                     |                          |              |
|--|-----------------------|-----------------------------|--------------------------|--------------|-----------------------------|--------------------------|--------------|-----------------------------|--------------------------|--------------|
|  |                       | Amt Found In Sample (ng/mL) | Amount Recovered (ng/mL) | Recovery (%) | Amt Found In Sample (ng/mL) | Amount Recovered (ng/mL) | Recovery (%) | Amt Found In Sample (ng/mL) | Amount Recovered (ng/mL) | Recovery (%) |
| Reagent Spike A<br>0.1 ng/mL (Dataset 012510A)                           | 0.100                 | ND                          | 0.0917                   | 92           | ND                          | 0.101                    | 101          | ND                          | 0.115                    | 115          |
| Reagent Spike B<br>1.0 ng/mL (Dataset 012510A)                           | 1.00                  | ND                          | 1.03                     | 103          | ND                          | 1.15                     | 115          | ND                          | 1.11                     | 111          |
| 25-1 Spike C<br>(L19874-1 Spk C, 1.0 ng/mL Lab Spike)<br>Dataset 012510A | 1.00                  | ND                          | 0.831                    | 83           | 0.181                       | 1.29                     | 111          | ND                          | 1.27                     | 127          |
| Reagent Spike A<br>0.1 ng/mL (Dataset 012810A)                           | 0.100                 | N/A                         | N/A                      | N/A          | N/A                         | N/A                      | N/A          | N/A                         | N/A                      | N/A          |
| Reagent Spike B<br>1.0 ng/mL (Dataset 012810A)                           | 1.00                  | N/A                         | N/A                      | N/A          | N/A                         | N/A                      | N/A          | N/A                         | N/A                      | N/A          |
| 25-1 Spike C<br>(L19784-1 Spk C, 1.0 ng/mL Lab Spike)<br>Dataset 012810A | 1.00                  | N/A                         | N/A                      | N/A          | N/A                         | N/A                      | N/A          | N/A                         | N/A                      | N/A          |

| Sample Description  | Amount Spiked (ng/mL) | C7 Acid                     |                          |              | C8 Acid                     |                          |              | C9 Acid                     |                          |              |
|---|-----------------------|-----------------------------|--------------------------|--------------|-----------------------------|--------------------------|--------------|-----------------------------|--------------------------|--------------|
|   |                       | Amt Found In Sample (ng/mL) | Amount Recovered (ng/mL) | Recovery (%) | Amt Found In Sample (ng/mL) | Amount Recovered (ng/mL) | Recovery (%) | Amt Found In Sample (ng/mL) | Amount Recovered (ng/mL) | Recovery (%) |
| Reagent Spike A<br>0.1 ng/mL (Dataset 012510A)                            | 0.100                 | ND                          | 0.110                    | 110          | ND                          | 0.0900                   | 90           | ND                          | 0.115                    | 115          |
| Reagent Spike B<br>1.0 ng/mL (Dataset 012510A)                            | 1.00                  | ND                          | 1.16                     | 116          | ND                          | 1.00                     | 100          | ND                          | 1.20                     | 120          |
| 25-1 Spike C<br>(L19874-1 Spk C, 1.0 ng/mL Lab Spike)<br>Dataset 012510A  | 1.00                  | ND                          | 1.19                     | 119          | ND                          | 1.11                     | 111          | N/A                         | N/A                      | N/A          |
| Reagent Spike A<br>0.1 ng/mL (Dataset 012810A)                            | 0.100                 | N/A                         | N/A                      | N/A          | N/A                         | N/A                      | N/A          | ND                          | 0.0730                   | 73           |
| Reagent Spike B<br>1.0 ng/mL (Dataset 012810A)                            | 1.00                  | N/A                         | N/A                      | N/A          | N/A                         | N/A                      | N/A          | ND                          | 0.824                    | 82           |
| 25-1 Spike C<br>(L19874-1 Spk C, 1.0 ng/mL Lab Spike)<br>Dataset 012810A  | 1.00                  | N/A                         | N/A                      | N/A          | N/A                         | N/A                      | N/A          | ND                          | 0.0972                   | 10^          |
| 34-1 Spike D<br>(L19874-28 Spk D, 1.0 ng/mL Lab Spike)<br>Dataset 012810A | 1.00                  | N/A                         | N/A                      | N/A          | N/A                         | N/A                      | N/A          | ND                          | 0.837                    | 84           |
| 34-1 Spike E<br>(L19784-28 Spk E, 1.0 ng/mL Lab Spike)<br>Dataset 012810A | 1.00                  | N/A                         | N/A                      | N/A          | N/A                         | N/A                      | N/A          | N/A                         | N/A                      | N/A          |

ND = Not detected = Response is below the LOD of 0.0125 ng/mL.

NQ = Not quantifiable = Response is between the LOD and the LOQ of 0.0250 ng/mL.

^ Confirmation analysis Laboratory Matrix Spike recovery is outside the acceptance criteria of 70 - 130% due to matrix effect

## Recovery Summary of Fluorochemical Residues in Water Samples

| Sample Description  | Amount Spiked (ng/mL) | C10 Acid                    |                          |              | C11 Acid                    |                          |              | C12 Acid                    |                          |              |
|---|-----------------------|-----------------------------|--------------------------|--------------|-----------------------------|--------------------------|--------------|-----------------------------|--------------------------|--------------|
|   |                       | Amt Found In Sample (ng/mL) | Amount Recovered (ng/mL) | Recovery (%) | Amt Found In Sample (ng/mL) | Amount Recovered (ng/mL) | Recovery (%) | Amt Found In Sample (ng/mL) | Amount Recovered (ng/mL) | Recovery (%) |
| Reagent Spike A<br>0.1 ng/mL (Dataset 012510A)                            | 0.100                 | ND                          | 0.108                    | 106          | ND                          | 0.109                    | 109          | ND                          | 0.117                    | 117          |
| Reagent Spike B<br>1.0 ng/mL (Dataset 012510A)                            | 1.00                  | ND                          | 1.05                     | 105          | ND                          | 1.23                     | 123          | ND                          | 1.30                     | 130          |
| 25-1 Spike C<br>(L19874-1 Spk C, 1.0 ng/mL Lab Spike)<br>Dataset 012510A  | 1.00                  | N/A                         | N/A                      | N/A          | N/A                         | N/A                      | N/A          | N/A                         | N/A                      | N/A          |
| Reagent Spike A<br>0.1 ng/mL (Dataset 012810A)                            | 0.100                 | ND                          | 0.0861                   | 86           | ND                          | 0.108                    | 108          | ND                          | 0.0963                   | 96           |
| Reagent Spike B<br>1.0 ng/mL (Dataset 012810A)                            | 1.00                  | ND                          | 0.922                    | 92           | ND                          | 1.00                     | 100          | ND                          | 1.03                     | 103          |
| 25-1 Spike C<br>(L19874-1 Spk C, 1.0 ng/mL Lab Spike)<br>Dataset 012810A  | 1.00                  | ND                          | 1.01                     | 101*         | ND                          | 2.52                     | 252*         | ND                          | 0.754                    | 75*          |
| 34-1 Spike D<br>(L19874-28 Spk D, 1.0 ng/mL Lab Spike)<br>Dataset 012810A | 1.00                  | ND                          | 0.867                    | 87           | ND                          | 1.23                     | 123          | ND                          | 1.160                    | 116          |
| 34-1 Spike E<br>(L19784-28 Spk E, 1.0 ng/mL Lab Spike)<br>Dataset 012810A | 1.00                  | N/A                         | N/A                      | N/A          | N/A                         | N/A                      | N/A          | N/A                         | N/A                      | N/A          |

| Sample Description  | Amount Spiked (ng/mL) | C13 Acid                    |                          |              | C14 Acid                    |                          |              |
|---|-----------------------|-----------------------------|--------------------------|--------------|-----------------------------|--------------------------|--------------|
|   |                       | Amt Found In Sample (ng/mL) | Amount Recovered (ng/mL) | Recovery (%) | Amt Found In Sample (ng/mL) | Amount Recovered (ng/mL) | Recovery (%) |
| Reagent Spike A<br>0.1 ng/mL (Dataset 012510A)                            | 0.100                 | ND                          | 0.0943                   | 94           | ND                          | 0.115                    | 115          |
| Reagent Spike B<br>1.0 ng/mL (Dataset 012510A)                            | 1.00                  | ND                          | 1.16                     | 116          | ND                          | 1.17                     | 117          |
| 25-1 Spike C<br>(L19874-1 Spk C, 1.0 ng/mL Lab Spike)<br>Dataset 012510A  | 1.00                  | N/A                         | N/A                      | N/A          | N/A                         | N/A                      | N/A          |
| Reagent Spike A<br>0.1 ng/mL (Dataset 012810A)                            | 0.100                 | ND                          | 0.0945                   | 95           | ND                          | 0.0883                   | 88           |
| Reagent Spike B<br>1.0 ng/mL (Dataset 012810A)                            | 1.00                  | ND                          | 0.982                    | 98           | ND                          | 1.07                     | 107          |
| 25-1 Spike C<br>(L19874-1 Spk C, 1.0 ng/mL Lab Spike)<br>Dataset 012810A  | 1.00                  | ND                          | 1.57                     | 157*         | ND                          | 2.78                     | 278*         |
| 34-1 Spike D<br>(L19874-28 Spk D, 1.0 ng/mL Lab Spike)<br>Dataset 012810A | 1.00                  | ND                          | 1.20                     | 120          | ND                          | 1.43                     | 143*         |
| 34-1 Spike E<br>(L19784-28 Spk E, 1.0 ng/mL Lab Spike)<br>Dataset 012810A | 1.00                  | N/A                         | N/A                      | N/A          | ND                          | 1.25                     | 125          |

ND = Not detected = Response is below the LOD of 0.0125 ng/mL.

\* Not quantifiable = Response is between the LOD and the LOQ of 0.0250 ng/mL.

Confirmation analysis Laboratory Matrix Spike recovery is outside the acceptance criteria of 70 - 130% due to matrix effect

\* Confirmation analysis IS recovery is <50% of the average standard IS value due to matrix effect. Value should be considered an estimate

## Recovery Summary of Fluorochemical Residues in Water Samples

| Sample Description  | Amount Spiked (ng/mL) | MeFOSAA                     |                          |              | EtFOSAA                     |                          |                  |
|---|-----------------------|-----------------------------|--------------------------|--------------|-----------------------------|--------------------------|------------------|
|   |                       | Amt Found in Sample (ng/mL) | Amount Recovered (ng/mL) | Recovery (%) | Amt Found in Sample (ng/mL) | Amount Recovered (ng/mL) | Recovery (%)     |
| Reagent Spike A<br>0.1 ng/mL (Dataset: 012510B)                           | 0.100                 | ND                          | 0.0984                   | 98           | ND                          | 0.123                    | 123              |
| Reagent Spike B<br>1.0 ng/mL (Dataset: 012510B)                           | 1.00                  | ND                          | 1.00                     | 100          | ND                          | 1.05                     | 105              |
| Reagent Spike A<br>0.1 ng/mL (Dataset: 012810B)                           | 0.100                 | N/A                         | N/A                      | N/A          | ND                          | 0.112                    | 112              |
| Reagent Spike B<br>1.0 ng/mL (Dataset: 012810B)                           | 1.00                  | N/A                         | N/A                      | N/A          | ND                          | 1.12                     | 112              |
| 25-1 Spike C<br>(L19874-1 Spk C, 0.1 ng/mL Lab Spike)<br>Dataset: 012510B | 0.100                 | ND                          | 0.123                    | 123          | N/A                         | N/A                      | N/A              |
| 25-1 Spike C<br>(L19874-1 Spk C, 1.0 ng/mL Lab Spike)<br>Dataset: 012810B | 1.00                  | N/A                         | N/A                      | N/A          | ND                          | 1.44                     | 144 <sup>^</sup> |

| Sample Description  | Amount Spiked (ng/mL) | MeFOSE                      |                          |              | EtFOSE                      |                          |              |
|---|-----------------------|-----------------------------|--------------------------|--------------|-----------------------------|--------------------------|--------------|
|   |                       | Amt Found in Sample (ng/mL) | Amount Recovered (ng/mL) | Recovery (%) | Amt Found in Sample (ng/mL) | Amount Recovered (ng/mL) | Recovery (%) |
| Reagent Spike A<br>0.1 ng/mL (Dataset: 012510B)                           | 0.100                 | ND                          | 0.109                    | 109          | ND                          | 0.118                    | 118          |
| Reagent Spike B<br>1.0 ng/mL (Dataset: 012510B)                           | 1.00                  | ND                          | 1.06                     | 106          | ND                          | 1.09                     | 109          |
| Reagent Spike A<br>0.1 ng/mL (Dataset: 012810B)                           | 0.100                 | ND                          | 0.110                    | 110          | ND                          | 0.109                    | 109          |
| Reagent Spike B<br>1.0 ng/mL (Dataset: 012810B)                           | 1.00                  | ND                          | 1.03                     | 103          | ND                          | 1.07                     | 107          |
| 25-1 Spike C<br>(L19874-1 Spk C, 0.1 ng/mL Lab Spike)<br>Dataset: 012510B | 0.100                 | N/A                         | N/A                      | N/A          | N/A                         | N/A                      | N/A          |
| 25-1 Spike C<br>(L19874-1 Spk C, 1.0 ng/mL Lab Spike)<br>Dataset: 012810B | 1.00                  | 2.60                        | 3.33                     | 73           | NQ                          | 0.695                    | 70           |

ND = Not detected = Response is below the LOD of 0.0125 ng/mL.

NQ = Not quantifiable = Response is between the LOD and the LOQ of 0.0250 ng/mL.

<sup>^</sup> Confirmation analysis of Laboratory Matrix Spike recovery is outside the acceptance criteria of 70 - 130% from matrix effect.



---

**Fluorotelomer Analysis by GC/MS**

**Sample ID: Method Blank**

**Date of Extraction: 1/12/2010**

**Date Analyzed: 1/13/2010**

| Analyte   | Amount Found (ng/mL) | LOQ (ng/mL) |
|-----------|----------------------|-------------|
| 7-2s FTOH | ND                   | 1.00        |
| 6-2 FTOH  | ND                   | 1.00        |
| 8-2 FTOH  | ND                   | 1.00        |
| 10-2 FTOH | ND                   | 1.00        |



3058 Research Drive  
State College, Pennsylvania 16801 USA  
Telephone: 814.272.1039  
Fax: 814.272.1019

## Analytical Report

---

### Fluorotelomer Analysis by GC/MS

**Sample ID: Lab Control Spike**

**Date of Extraction: 1/12/2010**

**Date Analyzed: 1/13/2010**

| Analyte   | Amount Found (ng/mL) | Amount Spiked (ng/mL) | LOQ (ng/mL) | % Recovery |
|-----------|----------------------|-----------------------|-------------|------------|
| 7-2s FTOH | 3.30                 | 5.00                  | 1.00        | 66         |
| 6-2 FTOH  | 8.20                 | 9.88                  | 1.00        | 83         |
| 8-2 FTOH  | 3.60                 | 5.11                  | 1.00        | 71         |
| 10-2 FTOH | 3.60                 | 5.03                  | 1.00        | 72         |





**MPI**  
RESEARCH

3058 Research Drive  
State College, Pennsylvania 16801 USA  
Telephone: 814.272.1039  
Fax: 814.272.1019

## Analytical Report

---

### Fluorotelomer Analysis by GC/MS

**Sample ID: Lab Control Spike Duplicate**

**Date of Extraction: 1/12/2010**

**Date Analyzed: 1/13/2010**

| Analyte   | Amount Found (ng/mL) | Amount Spiked (ng/mL) | LOQ (ng/mL) | % Recovery | % RPD |
|-----------|----------------------|-----------------------|-------------|------------|-------|
| 7-2s FTOH | 3.40                 | 5.00                  | 1.00        | 68         | 3.0   |
| 6-2 FTOH  | 8.50                 | 9.88                  | 1.00        | 86         | 3.6   |
| 8-2 FTOH  | 4.00                 | 5.11                  | 1.00        | 78         | 10.5  |
| 10-2 FTOH | 4.00                 | 5.03                  | 1.00        | 80         | 10.5  |

---

**Fluorotelomer Analysis by GC/MS**

**Sample ID: Method Blank**  
**Date of Extraction: 1/19/2010**  
**Date Analyzed: 1/19/2010**

| Analyte   | Amount Found (ng/mL) | LOQ (ng/mL) |
|-----------|----------------------|-------------|
| 7-2s FTOH | ND                   | 1.00        |
| 6-2 FTOH  | ND                   | 1.00        |
| 8-2 FTOH  | ND                   | 1.00        |
| 10-2 FTOH | ND                   | 1.00        |



**MPI**  
RESEARCH

3058 Research Drive  
State College, Pennsylvania 16801 USA  
Telephone: 814.272.1039  
Fax: 814.272.1019

## Analytical Report

---

### Fluorotelomer Analysis by GC/MS

**Sample ID: Lab Contol Spike**

**Date of Extraction: 1/19/2010**

**Date Analyzed: 1/19/2010**

| Analyte   | Amount Found (ng/mL) | Amount Spiked (ng/mL) | LOQ (ng/mL) | % Recovery |
|-----------|----------------------|-----------------------|-------------|------------|
| 7-2s FTOH | 3.40                 | 5.00                  | 1.00        | 68         |
| 6-2 FTOH  | 3.40                 | 4.94                  | 1.00        | 69         |
| 8-2 FTOH  | 3.30                 | 5.11                  | 1.00        | 65         |
| 10-2 FTOH | 3.40                 | 5.03                  | 1.00        | 68         |



3058 Research Drive  
State College, Pennsylvania 16801 USA  
Telephone: 814.272.1039  
Fax: 814.272.1019

## Analytical Report

### Fluorotelomer Analysis by GC/MS

**Sample ID: Lab Control Spike Duplicate**  
**Date of Extraction: 1/19/2010**  
**Date Analyzed: 1/19/2010**

| Analyte   | Amount Found (ng/mL) | Amount Spiked (ng/mL) | LOQ (ng/mL) | % Recovery | % RPD |
|-----------|----------------------|-----------------------|-------------|------------|-------|
| 7-2s FTOH | 3.70                 | 5.00                  | 1.00        | 74         | 8.5   |
| 6-2 FTOH  | 3.50                 | 4.94                  | 1.00        | 71         | 2.9   |
| 8-2 FTOH  | 3.30                 | 5.11                  | 1.00        | 65         | 0.0   |
| 10-2 FTOH | 3.60                 | 5.03                  | 1.00        | 72         | 5.7   |

## Recovery Summary of Fluorochemical Residues in Water Samples

| Sample Description  | Amount Spiked (ng/mL) | PFBS                        |                          |                  | PFHS                        |                          |                  | PFOS                        |                          |                 |
|---|-----------------------|-----------------------------|--------------------------|------------------|-----------------------------|--------------------------|------------------|-----------------------------|--------------------------|-----------------|
|   |                       | Amt Found in Sample (ng/mL) | Amount Recovered (ng/mL) | Recovery (%)     | Amt Found in Sample (ng/mL) | Amount Recovered (ng/mL) | Recovery (%)     | Amt Found in Sample (ng/mL) | Amount Recovered (ng/mL) | Recovery (%)    |
| Reagent Spike A<br>0.1 ng/mL (Dataset 012510A)                            | 0.100                 | ND                          | 0.0970                   | 97               | ND                          | 0.100                    | 100              | ND                          | 0.0994                   | 99              |
| Reagent Spike B<br>1.0 ng/mL (Dataset 012510A)                            | 1.00                  | ND                          | 1.09                     | 109              | ND                          | 1.08                     | 108              | ND                          | 1.10                     | 110             |
| 25-1 Spike C<br>(L19874-1 Spk C, 1.0 ng/mL Lab Spike)<br>Dataset 012510A  | 1.00                  | N/A                         | N/A                      | N/A              | N/A                         | N/A                      | N/A              | N/A                         | N/A                      | N/A             |
| Reagent Spike A<br>0.1 ng/mL (Dataset 012810A)                            | 0.100                 | ND                          | 0.0781                   | 78               | ND                          | 0.0731                   | 73               | ND                          | 0.0807                   | 81              |
| Reagent Spike B<br>1.0 ng/mL (Dataset 012810A)                            | 1.00                  | ND                          | 0.832                    | 83               | ND                          | 0.849                    | 85               | ND                          | 0.823                    | 82              |
| 25-1 Spike C<br>(L19874-1 Spk C, 1.0 ng/mL Lab Spike)<br>Dataset 012810A  | 1.00                  | ND                          | 8.04                     | 804 <sup>^</sup> | ND                          | 7.71                     | 771 <sup>^</sup> | ND                          | 0.889                    | 89 <sup>^</sup> |
| 34-1 Spike D<br>(L19874-28 Spk D, 1.0 ng/mL Lab Spike)<br>Dataset 012810A | 1.00                  | ND                          | 0.844                    | 84               | ND                          | 0.928                    | 93               | ND                          | 0.849                    | 85              |
| 34-1 Spike E<br>(L19784-28 Spk E, 1.0 ng/mL Lab Spike)<br>Dataset 012810A | 1.00                  | N/A                         | N/A                      | N/A              | N/A                         | N/A                      | N/A              | N/A                         | N/A                      | N/A             |

| Sample Description  | Amount Spiked (ng/mL) | FOSA                        |                          |                  | PFHpS                       |                          |                  | PFDS                        |                          |                  |
|---|-----------------------|-----------------------------|--------------------------|------------------|-----------------------------|--------------------------|------------------|-----------------------------|--------------------------|------------------|
|   |                       | Amt Found in Sample (ng/mL) | Amount Recovered (ng/mL) | Recovery (%)     | Amt Found in Sample (ng/mL) | Amount Recovered (ng/mL) | Recovery (%)     | Amt Found in Sample (ng/mL) | Amount Recovered (ng/mL) | Recovery (%)     |
| Reagent Spike A<br>0.1 ng/mL (Dataset 012510A)                            | 0.100                 | ND                          | 0.107                    | 107              | ND                          | 0.113                    | 113              | ND                          | 0.107                    | 107              |
| Reagent Spike B<br>1.0 ng/mL (Dataset 012510A)                            | 1.00                  | ND                          | 1.24                     | 124              | ND                          | 1.19                     | 119              | ND                          | 1.20                     | 120              |
| 25-1 Spike C<br>(L19874-1 Spk C, 1.0 ng/mL Lab Spike)<br>Dataset 012510A  | 1.00                  | N/A                         | N/A                      | N/A              | N/A                         | N/A                      | N/A              | N/A                         | N/A                      | N/A              |
| Reagent Spike A<br>0.1 ng/mL (Dataset 012810A)                            | 0.100                 | ND                          | 0.0754                   | 75               | ND                          | 0.0764                   | 76               | ND                          | 0.0806                   | 81               |
| Reagent Spike B<br>1.0 ng/mL (Dataset 012810A)                            | 1.00                  | ND                          | 0.927                    | 93               | ND                          | 0.859                    | 86               | ND                          | 0.869                    | 87               |
| 25-1 Spike C<br>(L19874-1 Spk C, 1.0 ng/mL Lab Spike)<br>Dataset 012810A  | 1.00                  | ND                          | 6.67                     | 667 <sup>^</sup> | ND                          | 2.64                     | 264 <sup>^</sup> | ND                          | 5.71                     | 571 <sup>^</sup> |
| 34-1 Spike D<br>(L19874-28 Spk D, 1.0 ng/mL Lab Spike)<br>Dataset 012810A | 1.00                  | ND                          | 0.957                    | 96               | ND                          | 0.887                    | 89               | ND                          | 0.862                    | 86               |
| 34-1 Spike E<br>(L19784-28 Spk E, 1.0 ng/mL Lab Spike)<br>Dataset 012810A | 1.00                  | N/A                         | N/A                      | N/A              | N/A                         | N/A                      | N/A              | N/A                         | N/A                      | N/A              |

ND = Not detected = Response is below the LOD of 0.0125 ng/mL.

NQ = Not quantifiable = Response is between the LOD and the LOQ of 0.0250 ng/mL.

Confirmation analysis Laboratory Matrix Spike recovery is outside the acceptance criteria of 70 - 130% due to matrix effect  
Confirmation analysis IS recovery is <50% of the average standard IS value due to matrix effect. Value should be considered an estimate

## Fluorotelomer Analysis by GC/MS

**Sample ID: L0019874-33 Matrix Spike; 35-3**

**Date of Extraction: 1/19/2010**

**Date Analyzed: 1/19/2010**

| Analyte   | Amount Found (ng/mL) | Amount Spiked (ng/mL) | LOQ (ng/mL) | % Recovery |
|-----------|----------------------|-----------------------|-------------|------------|
| 7-2s FTOH | 2.90                 | 5.00                  | 1.00        | 58         |
| 6-2 FTOH  | 2.70                 | 4.94                  | 1.00        | 55         |
| 8-2 FTOH  | 2.70                 | 5.11                  | 1.00        | 53         |
| 10-2 FTOH | 2.70                 | 5.03                  | 1.00        | 54         |

## Recovery Summary of <sup>13</sup>C PFOA (m+4) in Water Samples

| Client<br>Sample ID | MPI<br>Sample ID | Amount<br>Spiked<br>(ng/mL, ppb) | Amount<br>Recovered<br>(ng/mL, ppb) | Recovery<br>(%) |
|---------------------|------------------|----------------------------------|-------------------------------------|-----------------|
| N/A                 | Reagent Control  | 1.00                             | 1.01                                | 101             |
| N/A                 | Reagent Spike A  | 0.100                            | 0.0834                              | 83              |
| N/A                 | Reagent Spike B  | 1.00                             | 0.972                               | 97              |
| 25-1 Spike C        | L19874-1 Spk C   | 1.00                             | 1.08                                | 108             |
| 25-1                | L19874-1         | 1.00                             | 0.979                               | 98              |
| 26-1                | L19874-4         | 1.00                             | 0.962                               | 96              |
| 27-1                | L19874-7         | 1.00                             | 1.05                                | 105             |
| 28-1                | L19874-10        | 1.00                             | 1.17                                | 117             |
| 29-1                | L19874-13        | 1.00                             | 1.07                                | 107             |
| 30-1                | L19874-16        | 1.00                             | 1.09                                | 109             |
| 31-1                | L19874-19        | 1.00                             | 1.08                                | 108             |
| 32-1                | L19874-22        | 1.00                             | 1.08                                | 108             |
| 33-1                | L19874-25        | 1.00                             | 1.12                                | 112             |
| 34-1                | L19874-28        | 1.00                             | 1.09                                | 109             |
| 35-1                | L19874-31        | 1.00                             | 1.05                                | 105             |
| 4D1 Duplicate       | L19874-34        | 1.00                             | 1.05                                | 105             |

## Recovery Summary of 8:1 FTOH in Water Samples

| Client<br>Sample ID | MPI<br>Sample ID       | Amount<br>Spiked<br>(ng/mL, ppb) | Amount<br>Recovered<br>(ng/mL, ppb) | Recovery<br>(%) |
|---------------------|------------------------|----------------------------------|-------------------------------------|-----------------|
| N/A                 | Method Blank (1/12/10) | 5.00                             | 3.60                                | 72              |
| N/A                 | LCS (1/12/10)          | 5.00                             | 3.30                                | 66              |
| N/A                 | LCSD (1/12/10)         | 5.00                             | 3.70                                | 74              |
| N/A                 | Method Blank (1/19/10) | 5.00                             | 3.30                                | 66              |
| N/A                 | LCS (1/19/10)          | 5.00                             | 2.90                                | 58              |
| N/A                 | LCSD (1/19/10)         | 5.00                             | 3.40                                | 68              |
| 25-3                | L19874-3               | 5.00                             | 3.50                                | 70              |
| 26-3                | L19847-6               | 5.00                             | 2.50                                | 50              |
| 27-3                | L19874-9               | 5.00                             | 3.00                                | 60              |
| 28-3                | L19874-12              | 5.00                             | 2.10                                | 42              |
| 29-3                | L19874-15              | 5.00                             | 3.60                                | 72              |
| 30-3                | L19874-18              | 5.00                             | 2.70                                | 54              |
| 31-3                | L19874-21              | 5.00                             | 2.50                                | 50              |
| 32-3                | L19874-24              | 5.00                             | 0.80                                | 16              |
| 33-3                | L19874-27              | 5.00                             | 3.60                                | 72              |
| 34-3                | L19874-30              | 5.00                             | 3.80                                | 76              |
| 35-3                | L19874-33              | 5.00                             | 2.40                                | 48              |
| 35-3                | L19874-33 Duplicate    | 5.00                             | 2.20                                | 44              |
| 35-3                | L19874-33 Matrix Spike | 5.00                             | 2.60                                | 52              |
| 4D3 Duplicate       | L19874-36              | 5.00                             | 3.10                                | 62              |