

## ANALYTICAL REPORT

Job Number: 280-105698-1

Job Description: FAY-2018 Residential Sampling

For:  
Chemours Company FC, LLC The  
c/o AECOM  
Sabre Building, Suite 300  
4051 Ogletown Road  
Newark, DE 19713  
Attention: Michael Aucoin



Approved for release.  
Carissa N Cumine  
Project Manager II  
2/2/2018 5:08 PM

---

Designee for  
Michelle A Johnston, Project Manager II  
4955 Yarrow Street, Arvada, CO, 80002  
(303)736-0110  
michelle.johnston@testamericainc.com  
02/02/2018  
Revision: 1

cc: Barbara McGraw  
Kelly Rinehimer

The test results in this report relate only to the samples in this report and meet all requirements of NELAC, with any exceptions noted. Pursuant to NELAP, this report shall not be reproduced except in full, without the written approval of the laboratory. All questions regarding this report should be directed to the TestAmerica Denver Project Manager.

The Lab Certification ID# is 4025.

Reporting limits are adjusted for sample size used, dilutions and moisture content if applicable.

**TestAmerica Laboratories, Inc.**

TestAmerica Denver 4955 Yarrow Street, Arvada, CO 80002  
Tel (303) 736-0100 Fax (303) 431-7171 [www.testamericainc.com](http://www.testamericainc.com)

# Table of Contents

Cover Title Page . . . . .	1
Data Summaries . . . . .	4
Definitions . . . . .	4
Case Narrative . . . . .	5
Detection Summary . . . . .	7
Client Sample Results . . . . .	19
Default Detection Limits . . . . .	29
Surrogate Summary . . . . .	30
QC Sample Results . . . . .	31
QC Association . . . . .	33
Chronicle . . . . .	34
Certification Summary . . . . .	37
Method Summary . . . . .	38
Sample Summary . . . . .	39
Manual Integration Summary . . . . .	40
Reagent Traceability . . . . .	42
COAs . . . . .	48
Organic Sample Data . . . . .	74
LCMS . . . . .	74
8321A_HFPO_Du . . . . .	74
8321A_HFPO_Du QC Summary . . . . .	75
8321A_HFPO_Du Sample Data . . . . .	82
Standards Data . . . . .	124
8321A_HFPO_Du ICAL Data . . . . .	124
8321A_HFPO_Du CCAL Data . . . . .	167
Raw QC Data . . . . .	188

# Table of Contents

8321A_HFPO_Du Tune Data .....	188
8321A_HFPO_Du Blank Data .....	193
8321A_HFPO_Du LCS/LCSD Data .....	201
8321A_HFPO_Du MS/MSD Data .....	218
8321A_HFPO_Du Duplicate/Triplicate Data .....	222
8321A_HFPO_Du Run Logs .....	226
8321A_HFPO_Du Prep Data .....	229
Shipping and Receiving Documents .....	233
Client Chain of Custody .....	234
Sample Receipt Checklist .....	237

## Definitions/Glossary

Client: Chemours Company FC, LLC The  
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 280-105698-1

### Glossary

Abbreviation	These commonly used abbreviations may or may not be present in this report.
α	Listed under the "D" column to designate that the result is reported on a dry weight basis
%R	Percent Recovery
CFL	Contains Free Liquid
CNF	Contains No Free Liquid
DER	Duplicate Error Ratio (normalized absolute difference)
Dil Fac	Dilution Factor
DL	Detection Limit (DoD/DOE)
DL, RA, RE, IN	Indicates a Dilution, Re-analysis, Re-extraction, or additional Initial metals/anion analysis of the sample
DLC	Decision Level Concentration (Radiochemistry)
EDL	Estimated Detection Limit (Dioxin)
LOD	Limit of Detection (DoD/DOE)
LOQ	Limit of Quantitation (DoD/DOE)
MDA	Minimum Detectable Activity (Radiochemistry)
MDC	Minimum Detectable Concentration (Radiochemistry)
MDL	Method Detection Limit
ML	Minimum Level (Dioxin)
NC	Not Calculated
ND	Not Detected at the reporting limit (or MDL or EDL if shown)
PQL	Practical Quantitation Limit
QC	Quality Control
RER	Relative Error Ratio (Radiochemistry)
RL	Reporting Limit or Requested Limit (Radiochemistry)
RPD	Relative Percent Difference, a measure of the relative difference between two points
TEF	Toxicity Equivalent Factor (Dioxin)
TEQ	Toxicity Equivalent Quotient (Dioxin)



**CASE NARRATIVE**  
**Client: The Chemours Company FC, LLC**  
**Project: FAY-2018 Residential Sampling**  
**Report Number: 280-105698-1**

With the exceptions noted as flags or footnotes, standard analytical protocols were followed in the analysis of the samples and no problems were encountered or anomalies observed. In addition all laboratory quality control samples were within established control limits, with any exceptions noted below. Each sample was analyzed to achieve the lowest possible reporting limit within the constraints of the method. In some cases, due to interference or analytes present at high concentrations, samples were diluted. For diluted samples, the reporting limits are adjusted relative to the dilution required.

Calculations are performed before rounding to avoid round-off errors in calculated results.

All holding times were met and proper preservation noted for the methods performed on these samples, unless otherwise detailed in the individual sections below.

For samples requiring analysis at a dilution, the dilution factor has been multiplied by the Method Detection Limit (MDL) for each analyte and evaluated versus the project-specific reporting limit (PSRL). If the obtained value is below the PSRL, then the PSRL is preserved as the reporting limit for the diluted result, otherwise, the obtained value becomes the reporting limit. This is done in order to maintain the PSRL to meet project requirements at the request of the client and to report the lowest possible RL for each analyte.

**Revision - 2/2/2018**

The sample ID for FAY-D-7145BVTLE-W1-1-012218 (280-105698-7) was revised to FAY-D-7145BUTLE-W1-1-012218.

**Receipt**

The samples were received on 1/23/2018 9:30 AM; the samples arrived in good condition, properly preserved and, where required, on ice. The temperature of the cooler at receipt was 0.9° C.

**Receipt Exceptions**

The sample ID on the container labels of sample volume received for FAY-D-7012NC87-W1-1-012218 (280-105698-10) did not match the information on the chain of custody (COC). The container labels list a sample ID of FAY-D-7012NC87H-W1-1-012218 while the COC lists a sample ID of FAY-D-7012NC87-W1-1-012218. The laboratory logged the sample ID per the chain of custody. The client was notified on 1/23/2018. The laboratory revised the sample ID to match the revised chain of custody submitted to the laboratory by the client via email on 1/26/2018. The revised chain of custody is included in the final report.

The custody seal on the cooler, as it was observed at the time of sample receipt, was placed across the top of the cooler and underneath the shipping tape instead of across the seal of the cooler to demonstrate the absence of tampering with sample volume during transit. However, it can be noted that the shipping tape was in tact and there was no observable evidence of tampering with the sample volume during transit. The client was notified on 1/23/2018.

No other anomalies were observed during sample receipt.

**Standards**

Analytical standards were prepared using the acid form of the compound Perfluoro(2-propoxypropanoic) acid (HFPO-DA).

The surrogate compound, 13C3 HFPO-DA was introduced at the extraction step and was used as an internal standard for quantitation of HFPO-DA. The concentration of the surrogate spike is 0.2ug/L in water samples or 50ug/kg in soil samples.

**Sample Extraction and Analysis**

The samples presented in this report were extracted for the target analyte by TestAmerica Denver's SOP DV-OP-0019, Rev. 8 and analyzed for the target analyte by TestAmerica Denver's SOP DV-LC-0012, Rev. 14, with the exceptions of the items indicated in the DuPont QAS. Sample FAY-D-6394CHKFT-W1-1-012218 (280-105698-1) was chosen to be analyzed as a duplicate and also to be spiked with the target analyte.

For water samples a 250mL aliquot of each sample is extracted using solid phase extraction technique with methanol conditioned Weak Anion Exchange cartridges. Each sample is spiked with the internal standard/surrogate, prior to extraction. After the sample is passed through the cartridge, the analytes are eluted with 2%Formic Acid, 6mLs of HPLC grade MeOH and then with 4mL of 10% ammonium hydroxide in methanol. The final volume is brought to 5mL using reagent water and the extract is analyzed by LC/MS/MS.

The target analyte is separated from other components on a high-performance liquid chromatography (HPLC) C18 column with a mobile phase mixture of water containing 0.1% ammonium acetate and methanol. The mass spectrometer detector is operated in the electrospray (ESI) negative ion mode. The instrument is calibrated at 7 concentration levels (0.2, 0.5, 1.0, 2.0, 5.0, 10 and 20ug/L). The target analyte is detected as the perfluoro(2-propoxypropanoic) acid with the parent ion of 328.8 amu. The daughter ions used for analysis by LC/MS/MS are at 284.8 amu. The ratio of the peak areas to the two ions must be  $\pm 20\%$  of the ion ratios in the mid-point ICAL for qualitative identification. Sample results are quantitated using the internal standard dilution.

**Tuning and Calibration**

The instrument is tuned with a solution of the target analyte such that mass assignments are within  $\pm 0.5$  amu of the daughter ions. The instrument is calibrated with seven concentration levels from 0.2ug/L to 20ug/L. Linear regression ( $y=ax+b$ ) or quadratic functions ( $y=ax^2+bx+c$ ) are used with a correlation coefficient or coefficient of determination  $\geq 0.990$ . Following initial calibration (ICAL), an initial calibration blank (ICB) is tested, which consists of methanol spiked with the surrogate. The result for the target analyte must be less than one half the reporting limit (RL) to proceed.

Next an initial calibration verification (ICV) standard is tested. This is a mid-level concentration standard from a different vendor from the ICAL standard. If a different vendor is not available then, a different lot number from the same vendor is used. The ICV must be within 80-120% of the true value.

The quantitation limit verification standard is a standard from the same source as the ICAL tested run at the RL level to determine accuracy near the detection limit. This recovery must be within 70-130%.

Continuing calibration verification (CCV) standards are tested every 10 injections and are from the same source as the ICAL and are at mid-level concentration. The recovery of the CCVs must be 70-130% or recalibration is necessary.

#### **Method QC Samples**

The Method Blank is processed reagent water spiked with internal standard and prepared with each batch of 20 samples of the same matrix. All samples in the batch are processed at the same time and with the same reagents. The method blank must be less than the LOD or associated batch samples must be re-extracted and reanalyzed.

Each batch is prepared with a low- and a mid-level concentration spike Laboratory Control Samples (LCS). The recoveries of these samples must be within 70-130% or associated batch samples must be re-extracted and reanalyzed. If the recovery is biased high and samples are non-detect, results can be reported without re-extraction.

#### **Calculations**

##### **Sample Result Calculation**

For internal standard quantitation,

HFPO-DA Response = Area of HFPO-DA \* 13C3 HFPO-DA concentration / area of 13C3 HFPO-DA

Concentration in waters, ug/L = (Cex Vt)/(Vo)

Where:

Cex = Concentration measured in sample extract from the target analyte response (ng/mL)

Vt = Volume of total extract (mL)

Vo = Volume of water extracted (mL)

##### **2. Percent Recovery Calculation**

Spike Recovery = (SSR-SR)/(SA)x100%

Where:

SSR = Spike sample result

SR = Sample result

SA = Spike added

##### **3. Relative Percent Difference Calculation**

RPD = (SR - DR)/(1/2(SR+DR))x100

Where:

SR = Sample result

DR = Duplicate result

#### **HFPO-DA Analysis Anomalies**

Samples FAY-D-6394CHKFT-W1-1-012218 (280-105698-1), FAY-D-6394CHKFT-W1-1-012218-D (280-105698-2), FAY-D-6246CHKFT-W1-1-012218 (280-105698-3), FAY-D-318BOONE-W1-1-012218 (280-105698-4), FAY-D-41BOONE-W1-1-012218 (280-105698-5), FAY-D-FB-012218 (280-105698-6), FAY-D-7145BUTLE-W1-1-012218 (280-105698-7), FAY-D-1515SCLLY-W1-1-012218 (280-105698-8), FAY-D-7396SALIE-W1-1-012218 (280-105698-9) and FAY-D-7012NC87H-W1-1-012218 (280-105698-10) were analyzed for Perfluorinated Hydrocarbons in accordance with DV-LC-0012. The samples were prepared on 01/24/2018 and analyzed on 01/25/2018.

Calibration 9 (STD125) has been included in the raw data, but was not used in the Initial Calibration (ICAL).

Reporting limits have been adjusted accordingly for the initial volumes extracted.

No other analytical or quality issues were noted, other than those described in the Definitions/Glossary page.

## Fluorochemical Characterization of Water Samples Analytical Results

Chemours Sample Identification	TestAmerica Sample Identification	Collection Date/Time	Date Sample Received by TestAmerica	Analysis Date	HFPO-DA# (ug/L**)
FAY-D- 6394CHKFT-W1- 1-012218	280-105698-1	1/22/2018 16:19	1/23/2018	1/25/2018	0.033

# HFPO-DA – hexafluoropropylene oxide dimer acid, analyzed by method DV-LC-0012, Revision 14.

< = less than the stated value

\*\* ug/L – micrograms/liter (parts per billion)

### DEFINITIONS:

Reporting Limit (RL) for the procedure is approximately 0.010ug/L.

### RESULTS ARE CALCULATED ACCORDING TO THE FOLLOWING CRITERIA:

For samples analyzed in duplicate:

If the sample and laboratory duplicate are greater than 5X RL, the relative percent difference (RPD) is less than 20, the average value is reported. If the RPD is greater than 20, the higher value is reported.


If the sample or laboratory duplicate are less than 5X RL, and the absolute difference is less than RL, the average value is reported. If the absolute difference is greater than the RL, the higher value is reported.

Matrix Spike Recoveries:

Acceptable Range: 70%-130%

TestAmerica Sample ID	Matrix Spike Recoveries
280-105698-1	94%

### SUBMITTED BY:

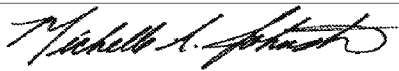

1/26/2018

Michelle A. Johnston, Project Manager

Date

## Fluorochemical Characterization of Water Samples Analytical Results

Chemours Sample Identification	TestAmerica Sample Identification	Collection Date/Time	Date Sample Received by TestAmerica	Analysis Date	HFPO-DA# (ug/L**)
FAY-D- 6394CHKFT-W1- 1-012218-D	280-105698-2	1/22/2018 16:19	1/23/2018	1/25/2018	0.032

# HFPO-DA – hexafluoropropylene oxide dimer acid, analyzed by method DV-LC-0012, Revision 14.

< = less than the stated value

\*\* ug/L – micrograms/liter (parts per billion)

### DEFINITIONS:

Reporting Limit (RL) for the procedure is approximately 0.010ug/L.

### RESULTS ARE CALCULATED ACCORDING TO THE FOLLOWING CRITERIA:

For samples analyzed in duplicate:

If the sample and laboratory duplicate are greater than 5X RL, the relative percent difference (RPD) is less than 20, the average value is reported. If the RPD is greater than 20, the higher value is reported.


If the sample or laboratory duplicate are less than 5X RL, and the absolute difference is less than RL, the average value is reported. If the absolute difference is greater than the RL, the higher value is reported.

Matrix Spike Recoveries:

Acceptable Range: 70%-130%

TestAmerica Sample ID	Matrix Spike Recoveries
280-105698-1	94%

### SUBMITTED BY:

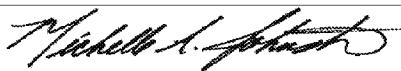

1/26/2018

Michelle A. Johnston, Project Manager

Date

## Fluorochemical Characterization of Water Samples Analytical Results

Chemours Sample Identification	TestAmerica Sample Identification	Collection Date/Time	Date Sample Received by TestAmerica	Analysis Date	HFPO-DA# (ug/L**)
FAY-D- 6246CHKFT-W1- 1-012218	280-105698-3	1/22/2018 14:56	1/23/2018	1/25/2018	0.052

# HFPO-DA – hexafluoropropylene oxide dimer acid, analyzed by method DV-LC-0012, Revision 14.

< = less than the stated value

\*\* ug/L – micrograms/liter (parts per billion)

### DEFINITIONS:

Reporting Limit (RL) for the procedure is approximately 0.010ug/L.

### RESULTS ARE CALCULATED ACCORDING TO THE FOLLOWING CRITERIA:

For samples analyzed in duplicate:

If the sample and laboratory duplicate are greater than 5X RL, the relative percent difference (RPD) is less than 20, the average value is reported. If the RPD is greater than 20, the higher value is reported.


If the sample or laboratory duplicate are less than 5X RL, and the absolute difference is less than RL, the average value is reported. If the absolute difference is greater than the RL, the higher value is reported.

Matrix Spike Recoveries:

Acceptable Range: 70%-130%

TestAmerica Sample ID	Matrix Spike Recoveries
280-105698-1	94%

### SUBMITTED BY:

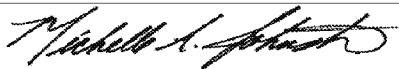

1/26/2018

Michelle A. Johnston, Project Manager

Date

## Fluorochemical Characterization of Water Samples Analytical Results

Chemours Sample Identification	TestAmerica Sample Identification	Collection Date/Time	Date Sample Received by TestAmerica	Analysis Date	HFPO-DA# (ug/L**)
FAY-D- 318BOONE-W1- 1-012218	280-105698-4	1/22/2018 14:15	1/23/2018	1/25/2018	0.044

# HFPO-DA – hexafluoropropylene oxide dimer acid, analyzed by method DV-LC-0012, Revision 14.

< = less than the stated value

\*\* ug/L – micrograms/liter (parts per billion)

### DEFINITIONS:

Reporting Limit (RL) for the procedure is approximately 0.010ug/L.

### RESULTS ARE CALCULATED ACCORDING TO THE FOLLOWING CRITERIA:

For samples analyzed in duplicate:

If the sample and laboratory duplicate are greater than 5X RL, the relative percent difference (RPD) is less than 20, the average value is reported. If the RPD is greater than 20, the higher value is reported.


If the sample or laboratory duplicate are less than 5X RL, and the absolute difference is less than RL, the average value is reported. If the absolute difference is greater than the RL, the higher value is reported.

Matrix Spike Recoveries:

Acceptable Range: 70%-130%

TestAmerica Sample ID	Matrix Spike Recoveries
280-105698-1	94%

### SUBMITTED BY:

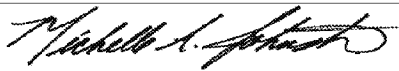

1/26/2018

Michelle A. Johnston, Project Manager

Date

## Fluorochemical Characterization of Water Samples Analytical Results

Chemours Sample Identification	TestAmerica Sample Identification	Collection Date/Time	Date Sample Received by TestAmerica	Analysis Date	HFPO-DA# (ug/L**)
FAY-D- 41BOONE-W1-1- 012218	280-105698-5	1/22/2018 13:49	1/23/2018	1/25/2018	<0.010

# HFPO-DA – hexafluoropropylene oxide dimer acid, analyzed by method DV-LC-0012, Revision 14.

< = less than the stated value

\*\* ug/L – micrograms/liter (parts per billion)

### DEFINITIONS:

Reporting Limit (RL) for the procedure is approximately 0.010ug/L.

### RESULTS ARE CALCULATED ACCORDING TO THE FOLLOWING CRITERIA:

For samples analyzed in duplicate:

If the sample and laboratory duplicate are greater than 5X RL, the relative percent difference (RPD) is less than 20, the average value is reported. If the RPD is greater than 20, the higher value is reported.


If the sample or laboratory duplicate are less than 5X RL, and the absolute difference is less than RL, the average value is reported. If the absolute difference is greater than the RL, the higher value is reported.

### Matrix Spike Recoveries:

Acceptable Range: 70%-130%

TestAmerica Sample ID	Matrix Spike Recoveries
280-105698-1	94%

### SUBMITTED BY:

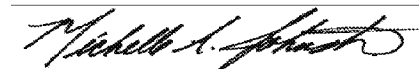

1/26/2018

Michelle A. Johnston, Project Manager

Date

## Fluorochemical Characterization of Water Samples Analytical Results

Chemours Sample Identification	TestAmerica Sample Identification	Collection Date/Time	Date Sample Received by TestAmerica	Analysis Date	HFPO-DA# (ug/L**)
FAY-D-FB- 012218	280-105698-6	1/22/2018 7:30	1/23/2018	1/25/2018	<0.010

# HFPO-DA – hexafluoropropylene oxide dimer acid, analyzed by method DV-LC-0012, Revision 14.

< = less than the stated value

\*\* ug/L – micrograms/liter (parts per billion)

### DEFINITIONS:

Reporting Limit (RL) for the procedure is approximately 0.010ug/L.

### RESULTS ARE CALCULATED ACCORDING TO THE FOLLOWING CRITERIA:

For samples analyzed in duplicate:

If the sample and laboratory duplicate are greater than 5X RL, the relative percent difference (RPD) is less than 20, the average value is reported. If the RPD is greater than 20, the higher value is reported.


If the sample or laboratory duplicate are less than 5X RL, and the absolute difference is less than RL, the average value is reported. If the absolute difference is greater than the RL, the higher value is reported.

### Matrix Spike Recoveries:

Acceptable Range: 70%-130%

TestAmerica Sample ID	Matrix Spike Recoveries
280-105698-1	94%

### SUBMITTED BY:

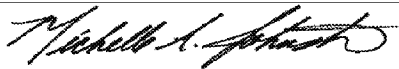

1/26/2018

Michelle A. Johnston, Project Manager

Date



## Fluorochemical Characterization of Water Samples Analytical Results

Chemours Sample Identification	TestAmerica Sample Identification	Collection Date/Time	Date Sample Received by TestAmerica	Analysis Date	HFPO-DA# (ug/L**)
FAY-D- 7145BUTLE-W1- 1-012218	280-105698-7	1/22/2018 14:03	1/23/2018	1/25/2018	0.080

# HFPO-DA – hexafluoropropylene oxide dimer acid, analyzed by method DV-LC-0012, Revision 14.

< = less than the stated value

\*\* ug/L – micrograms/liter (parts per billion)

### DEFINITIONS:

Reporting Limit (RL) for the procedure is approximately 0.010ug/L.

### RESULTS ARE CALCULATED ACCORDING TO THE FOLLOWING CRITERIA:

For samples analyzed in duplicate:

If the sample and laboratory duplicate are greater than 5X RL, the relative percent difference (RPD) is less than 20, the average value is reported. If the RPD is greater than 20, the higher value is reported.

If the sample or laboratory duplicate are less than 5X RL, and the absolute difference is less than RL, the average value is reported. If the absolute difference is greater than the RL, the higher value is reported.

### Matrix Spike Recoveries:

Acceptable Range: 70%-130%

TestAmerica Sample ID	Matrix Spike Recoveries
280-105698-1	94%

### SUBMITTED BY:

*Carissalumine*

2/2/2018

Carissa Cumine, Project Manager

Date

## Fluorochemical Characterization of Water Samples Analytical Results

Chemours Sample Identification	TestAmerica Sample Identification	Collection Date/Time	Date Sample Received by TestAmerica	Analysis Date	HFPO-DA# (ug/L**)
FAY-D- 1515SCLLY-W1- 1-012218	280-105698-8	1/22/2018 14:23	1/23/2018	1/25/2018	<0.010

# HFPO-DA – hexafluoropropylene oxide dimer acid, analyzed by method DV-LC-0012, Revision 14.

< = less than the stated value

\*\* ug/L – micrograms/liter (parts per billion)

### DEFINITIONS:

Reporting Limit (RL) for the procedure is approximately 0.010ug/L.

### RESULTS ARE CALCULATED ACCORDING TO THE FOLLOWING CRITERIA:

For samples analyzed in duplicate:

If the sample and laboratory duplicate are greater than 5X RL, the relative percent difference (RPD) is less than 20, the average value is reported. If the RPD is greater than 20, the higher value is reported.


If the sample or laboratory duplicate are less than 5X RL, and the absolute difference is less than RL, the average value is reported. If the absolute difference is greater than the RL, the higher value is reported.

Matrix Spike Recoveries:

Acceptable Range: 70%-130%

TestAmerica Sample ID	Matrix Spike Recoveries
280-105698-1	94%

### SUBMITTED BY:

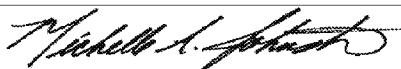

1/26/2018

Michelle A. Johnston, Project Manager

Date

## Fluorochemical Characterization of Water Samples Analytical Results

Chemours Sample Identification	TestAmerica Sample Identification	Collection Date/Time	Date Sample Received by TestAmerica	Analysis Date	HFPO-DA# (ug/L**)
FAY-D- 7396SALIE-W1- 1-012218	280-105698-9	1/22/2018 15:39	1/23/2018	1/25/2018	<0.010

# HFPO-DA – hexafluoropropylene oxide dimer acid, analyzed by method DV-LC-0012, Revision 14.

< = less than the stated value

\*\* ug/L – micrograms/liter (parts per billion)

### DEFINITIONS:

Reporting Limit (RL) for the procedure is approximately 0.010ug/L.

### RESULTS ARE CALCULATED ACCORDING TO THE FOLLOWING CRITERIA:

For samples analyzed in duplicate:

If the sample and laboratory duplicate are greater than 5X RL, the relative percent difference (RPD) is less than 20, the average value is reported. If the RPD is greater than 20, the higher value is reported.


If the sample or laboratory duplicate are less than 5X RL, and the absolute difference is less than RL, the average value is reported. If the absolute difference is greater than the RL, the higher value is reported.

Matrix Spike Recoveries:

Acceptable Range: 70%-130%

TestAmerica Sample ID	Matrix Spike Recoveries
280-105698-1	94%

### SUBMITTED BY:

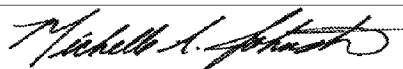

1/26/2018

Michelle A. Johnston, Project Manager

Date

## Fluorochemical Characterization of Water Samples Analytical Results

Chemours Sample Identification	TestAmerica Sample Identification	Collection Date/Time	Date Sample Received by TestAmerica	Analysis Date	HFPO-DA# (ug/L**)
FAY-D- 7012NC87H-W1- 1-012218	280-105698-10	1/22/2018 16:21	1/23/2018	1/25/2018	0.027

# HFPO-DA – hexafluoropropylene oxide dimer acid, analyzed by method DV-LC-0012, Revision 14.

< = less than the stated value

\*\* ug/L – micrograms/liter (parts per billion)

### DEFINITIONS:

Reporting Limit (RL) for the procedure is approximately 0.010ug/L.

### RESULTS ARE CALCULATED ACCORDING TO THE FOLLOWING CRITERIA:

For samples analyzed in duplicate:

If the sample and laboratory duplicate are greater than 5X RL, the relative percent difference (RPD) is less than 20, the average value is reported. If the RPD is greater than 20, the higher value is reported.


If the sample or laboratory duplicate are less than 5X RL, and the absolute difference is less than RL, the average value is reported. If the absolute difference is greater than the RL, the higher value is reported.

Matrix Spike Recoveries:

Acceptable Range: 70%-130%

TestAmerica Sample ID	Matrix Spike Recoveries
280-105698-1	94%

### SUBMITTED BY:

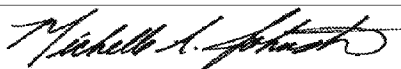

1/29/2018

Michelle A. Johnston, Project Manager

Date

## Executive Summary

Client: Chemours Company FC, LLC The

Job Number: 280-105698-1

### 8321A : HFPO-DA

Lab Sample ID	Client Sample ID	Analyte	Individual Result (ug/L)	Final Result (ug/L)	RL
280-105698-1	FAY-D-6394CHKFT-W1-1-012218	HFPO-DA	0.033	0.033	0.010
280-105698-1 DU	FAY-D-6394CHKFT-W1-1-012218	HFPO-DA	0.033		0.010
280-105698-2	FAY-D-6394CHKFT-W1-1-012218-D	HFPO-DA	0.032	0.032	0.010
280-105698-3	FAY-D-6246CHKFT-W1-1-012218	HFPO-DA	0.052	0.052	0.010
280-105698-4	FAY-D-318BOONE-W1-1-012218	HFPO-DA	0.044	0.044	0.010
280-105698-5	FAY-D-41BOONE-W1-1-012218	HFPO-DA	<0.010	<0.010	0.010
280-105698-6	FAY-D-FB-012218	HFPO-DA	<0.010	<0.010	0.010
280-105698-7	FAY-D-7145BUTLE-W1-1-012218	HFPO-DA	0.080	0.080	0.010
280-105698-8	FAY-D-1515SCLLY-W1-1-012218	HFPO-DA	<0.010	<0.010	0.010
280-105698-9	FAY-D-7396SALIE-W1-1-012218	HFPO-DA	<0.010	<0.010	0.010
280-105698-10	FAY-D-7012NC87H-W1-1-012218	HFPO-DA	0.027	0.027	0.010

(a) Method 8321A

(b) DUP or REP indicates a laboratory duplicate.

(c) If the sample and laboratory duplicate are both greater than 5X the RL and the relative percent difference (RPD) is less than 20, the average value is reported. If the RPD is greater than 20, the higher of the sample and laboratory duplicate value is reported. If the sample and/or laboratory duplicate are less than 5X the RL, and the absolute difference between the sample and laboratory duplicate is less than the RL, the average value is reported. If the absolute difference is greater than the RL, the higher of the sample and laboratory duplicate value is reported. If either the sample or the duplicate result is greater than or equal to the RL and the other is less than the RL, then the higher of the two is reported.

(d) Moisture Determined by ASTM D2216.

(e) Reporting Limit (RL) = The concentration equivalent to the low calibration standard.

## Detection Summary

Client: Chemours Company FC, LLC The  
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 280-105698-1

Client Sample ID: FAY-D-6394CHKFT-W1-1-012218

Lab Sample ID: 280-105698-1

Analyte	Result	Qualifier	RL	MDL	Unit	Dil Fac	D	Method	Prep Type
HFPO-DA	0.033		0.010		ug/L	1		8321A	Total/NA

Client Sample ID: FAY-D-6394CHKFT-W1-1-012218-D

Lab Sample ID: 280-105698-2

Analyte	Result	Qualifier	RL	MDL	Unit	Dil Fac	D	Method	Prep Type
HFPO-DA	0.032		0.010		ug/L	1		8321A	Total/NA

Client Sample ID: FAY-D-6246CHKFT-W1-1-012218

Lab Sample ID: 280-105698-3

Analyte	Result	Qualifier	RL	MDL	Unit	Dil Fac	D	Method	Prep Type
HFPO-DA	0.052		0.010		ug/L	1		8321A	Total/NA

Client Sample ID: FAY-D-318BOONE-W1-1-012218

Lab Sample ID: 280-105698-4

Analyte	Result	Qualifier	RL	MDL	Unit	Dil Fac	D	Method	Prep Type
HFPO-DA	0.044		0.010		ug/L	1		8321A	Total/NA

Client Sample ID: FAY-D-41BOONE-W1-1-012218

Lab Sample ID: 280-105698-5

No Detections.

Client Sample ID: FAY-D-FB-012218

Lab Sample ID: 280-105698-6

No Detections.

Client Sample ID: FAY-D-7145BUTLE-W1-1-012218

Lab Sample ID: 280-105698-7

Analyte	Result	Qualifier	RL	MDL	Unit	Dil Fac	D	Method	Prep Type
HFPO-DA	0.080		0.010		ug/L	1		8321A	Total/NA

Client Sample ID: FAY-D-1515SCLLY-W1-1-012218

Lab Sample ID: 280-105698-8

No Detections.

Client Sample ID: FAY-D-7396SALIE-W1-1-012218

Lab Sample ID: 280-105698-9

No Detections.

Client Sample ID: FAY-D-7012NC87H-W1-1-012218

Lab Sample ID: 280-105698-10

Analyte	Result	Qualifier	RL	MDL	Unit	Dil Fac	D	Method	Prep Type
HFPO-DA	0.027		0.010		ug/L	1		8321A	Total/NA

This Detection Summary does not include radiochemical test results.

TestAmerica Denver

# Client Sample Results

Client: Chemours Company FC, LLC The  
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 280-105698-1

Client Sample ID: FAY-D-6394CHKFT-W1-1-012218

Lab Sample ID: 280-105698-1

Date Collected: 01/22/18 16:19

Matrix: Water

Date Received: 01/23/18 09:30

## Method: 8321A - HFPO-DA

Analyte	Result	Qualifier	RL	MDL	Unit	D	Prepared	Analyzed	Dil Fac
HFPO-DA	0.033		0.010		ug/L		01/24/18 15:20	01/25/18 10:54	1
Surrogate	%Recovery	Qualifier	Limits				Prepared	Analyzed	Dil Fac
13C3 HFPO-DA	89		50 - 200				01/24/18 15:20	01/25/18 10:54	1

# Client Sample Results

Client: Chemours Company FC, LLC The  
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 280-105698-1

Client Sample ID: FAY-D-6394CHKFT-W1-1-012218-D

Lab Sample ID: 280-105698-2

Date Collected: 01/22/18 16:19

Matrix: Water

Date Received: 01/23/18 09:30

## Method: 8321A - HFPO-DA

Analyte	Result	Qualifier	RL	MDL	Unit	D	Prepared	Analyzed	Dil Fac
HFPO-DA	0.032		0.010		ug/L		01/24/18 15:20	01/25/18 11:07	1
Surrogate	%Recovery	Qualifier	Limits				Prepared	Analyzed	Dil Fac
13C3 HFPO-DA	88		50 - 200				01/24/18 15:20	01/25/18 11:07	1

TestAmerica Denver



# Client Sample Results

Client: Chemours Company FC, LLC The  
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 280-105698-1

Client Sample ID: FAY-D-6246CHKFT-W1-1-012218

Lab Sample ID: 280-105698-3

Date Collected: 01/22/18 14:56

Matrix: Water

Date Received: 01/23/18 09:30

## Method: 8321A - HFPO-DA

Analyte	Result	Qualifier	RL	MDL	Unit	D	Prepared	Analyzed	Dil Fac
HFPO-DA	0.052		0.010		ug/L		01/24/18 15:20	01/25/18 11:10	1
Surrogate	%Recovery	Qualifier	Limits				Prepared	Analyzed	Dil Fac
13C3 HFPO-DA	89		50 - 200				01/24/18 15:20	01/25/18 11:10	1

TestAmerica Denver

# Client Sample Results

Client: Chemours Company FC, LLC The  
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 280-105698-1

Client Sample ID: FAY-D-318BOONE-W1-1-012218

Lab Sample ID: 280-105698-4

Date Collected: 01/22/18 14:15

Matrix: Water

Date Received: 01/23/18 09:30

## Method: 8321A - HFPO-DA

Analyte	Result	Qualifier	RL	MDL	Unit	D	Prepared	Analyzed	Dil Fac
HFPO-DA	0.044		0.010		ug/L		01/24/18 15:20	01/25/18 11:14	1
Surrogate	%Recovery	Qualifier	Limits				Prepared	Analyzed	Dil Fac
13C3 HFPO-DA	92		50 - 200				01/24/18 15:20	01/25/18 11:14	1

TestAmerica Denver

# Client Sample Results

Client: Chemours Company FC, LLC The  
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 280-105698-1

Client Sample ID: FAY-D-41BOONE-W1-1-012218

Lab Sample ID: 280-105698-5

Date Collected: 01/22/18 13:49

Matrix: Water

Date Received: 01/23/18 09:30

## Method: 8321A - HFPO-DA

Analyte	Result	Qualifier	RL	MDL	Unit	D	Prepared	Analyzed	Dil Fac
HFPO-DA	<0.010		0.010		ug/L		01/24/18 15:20	01/25/18 11:17	1

Surrogate	%Recovery	Qualifier	Limits	Prepared	Analyzed	Dil Fac
13C3 HFPO-DA	88		50 - 200	01/24/18 15:20	01/25/18 11:17	1

TestAmerica Denver

# Client Sample Results

Client: Chemours Company FC, LLC The  
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 280-105698-1

Client Sample ID: FAY-D-FB-012218

Lab Sample ID: 280-105698-6

Date Collected: 01/22/18 07:30

Matrix: Water

Date Received: 01/23/18 09:30

## Method: 8321A - HFPO-DA

Analyte	Result	Qualifier	RL	MDL	Unit	D	Prepared	Analyzed	Dil Fac
HFPO-DA	<0.010		0.010		ug/L		01/24/18 15:20	01/25/18 11:20	1

Surrogate	%Recovery	Qualifier	Limits	Prepared	Analyzed	Dil Fac
13C3 HFPO-DA	92		50 - 200	01/24/18 15:20	01/25/18 11:20	1

TestAmerica Denver

# Client Sample Results

Client: Chemours Company FC, LLC The  
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 280-105698-1

Client Sample ID: FAY-D-7145BUTLE-W1-1-012218

Lab Sample ID: 280-105698-7

Date Collected: 01/22/18 14:03

Matrix: Water

Date Received: 01/23/18 09:30

## Method: 8321A - HFPO-DA

Analyte	Result	Qualifier	RL	MDL	Unit	D	Prepared	Analyzed	Dil Fac
HFPO-DA	0.080		0.010		ug/L		01/24/18 15:20	01/25/18 11:23	1
Surrogate	%Recovery	Qualifier	Limits				Prepared	Analyzed	Dil Fac
13C3 HFPO-DA	84		50 - 200				01/24/18 15:20	01/25/18 11:23	1

TestAmerica Denver

# Client Sample Results

Client: Chemours Company FC, LLC The  
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 280-105698-1

Client Sample ID: FAY-D-1515SCLLY-W1-1-012218

Lab Sample ID: 280-105698-8

Date Collected: 01/22/18 14:23

Matrix: Water

Date Received: 01/23/18 09:30

## Method: 8321A - HFPO-DA

Analyte	Result	Qualifier	RL	MDL	Unit	D	Prepared	Analyzed	Dil Fac
HFPO-DA	<0.010		0.010		ug/L		01/24/18 15:20	01/25/18 11:27	1

Surrogate	%Recovery	Qualifier	Limits	Prepared	Analyzed	Dil Fac
13C3 HFPO-DA	88		50 - 200	01/24/18 15:20	01/25/18 11:27	1

TestAmerica Denver

# Client Sample Results

Client: Chemours Company FC, LLC The  
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 280-105698-1

Client Sample ID: FAY-D-7396SALIE-W1-1-012218

Lab Sample ID: 280-105698-9

Date Collected: 01/22/18 15:39

Matrix: Water

Date Received: 01/23/18 09:30

## Method: 8321A - HFPO-DA

Analyte	Result	Qualifier	RL	MDL	Unit	D	Prepared	Analyzed	Dil Fac
HFPO-DA	<0.010		0.010		ug/L		01/24/18 15:20	01/25/18 11:33	1

Surrogate	%Recovery	Qualifier	Limits	Prepared	Analyzed	Dil Fac
13C3 HFPO-DA	86		50 - 200	01/24/18 15:20	01/25/18 11:33	1

TestAmerica Denver

# Client Sample Results

Client: Chemours Company FC, LLC The  
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 280-105698-1

Client Sample ID: FAY-D-7012NC87H-W1-1-012218

Lab Sample ID: 280-105698-10

Date Collected: 01/22/18 16:21

Matrix: Water

Date Received: 01/23/18 09:30

## Method: 8321A - HFPO-DA

Analyte	Result	Qualifier	RL	MDL	Unit	D	Prepared	Analyzed	Dil Fac
HFPO-DA	0.027		0.010		ug/L		01/24/18 15:20	01/25/18 11:36	1
Surrogate	%Recovery	Qualifier	Limits				Prepared	Analyzed	Dil Fac
13C3 HFPO-DA	86		50 - 200				01/24/18 15:20	01/25/18 11:36	1

TestAmerica Denver



## Default Detection Limits

Client: Chemours Company FC, LLC The  
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 280-105698-1

Method: 8321A - HFPO-DA

Prep: 3535

Analyte	RL	MDL	Units	Method
HFPO-DA	0.010	0.0051	ug/L	8321A

# Surrogate Summary

Client: Chemours Company FC, LLC The  
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 280-105698-1

Method: 8321A - HFPO-DA

Matrix: Water

Prep Type: Total/NA

Percent Surrogate Recovery (Acceptance Limits)		
Lab Sample ID	Client Sample ID	HFPODA (50-200)
280-105698-1	FAY-D-6394CHKFT-W1-1-0122	89
280-105698-1 DU	FAY-D-6394CHKFT-W1-1-0122	85
280-105698-1 MS	FAY-D-6394CHKFT-W1-1-0122	86
280-105698-2	FAY-D-6394CHKFT-W1-1-0122	88
280-105698-3	FAY-D-6246CHKFT-W1-1-0122	89
280-105698-4	FAY-D-318BOONE-W1-1-01221	92
280-105698-5	FAY-D-41BOONE-W1-1-01221E	88
280-105698-6	FAY-D-FB-012218	92
280-105698-7	FAY-D-7145BUTLE-W1-1-0122	84
280-105698-8	FAY-D-1515SCLLY-W1-1-01221	88
280-105698-9	FAY-D-7396SALIE-W1-1-01221	86
280-105698-10	FAY-D-7012NC87H-W1-1-0122	86
DLCK 280-390728/12	Lab Control Sample	102
LCS 280-402648/2-A	Lab Control Sample	91
LCSD 280-402648/3-A	Lab Control Sample Dup	90
LLCS 280-402648/4-A	Lab Control Sample	94
MB 280-402648/1-A	Method Blank	93
<b>Surrogate Legend</b>		
HFPODA = 13C3 HFPO-DA		

TestAmerica Denver

# QC Sample Results

Client: Chemours Company FC, LLC The  
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 280-105698-1

Method: 8321A - HFPO-DA

Lab Sample ID: DLCK 280-390728/12

Matrix: Water

Analysis Batch: 390728

Client Sample ID: Lab Control Sample

Prep Type: Total/NA

Analyte		Spike Added	DLCK Result	DLCK Qualifier	Unit	D	%Rec	% Rec. Limits
HFPO-DA		0.250	<0.50		ug/L		78	70 - 130
Surrogate	DLCK %Recovery	DLCK Qualifier	Limits					
13C3 HFPO-DA	102		50 - 200					

Lab Sample ID: MB 280-402648/1-A

Matrix: Water

Analysis Batch: 402806

Client Sample ID: Method Blank

Prep Type: Total/NA

Prep Batch: 402648

Analyte	MB Result	MB Qualifier	RL	MDL	Unit	D	Prepared	Analyzed	Dil Fac
HFPO-DA	<0.010		0.010		ug/L		01/24/18 15:20	01/25/18 10:41	1
Surrogate	MB %Recovery	MB Qualifier	Limits				Prepared	Analyzed	Dil Fac
13C3 HFPO-DA	93		50 - 200				01/24/18 15:20	01/25/18 10:41	1

Lab Sample ID: LCS 280-402648/2-A

Matrix: Water

Analysis Batch: 402806

Client Sample ID: Lab Control Sample

Prep Type: Total/NA

Prep Batch: 402648

Analyte		Spike Added	LCS Result	LCS Qualifier	Unit	D	%Rec	% Rec. Limits
HFPO-DA		0.200	0.193		ug/L		96	70 - 130
Surrogate	LCS %Recovery	LCS Qualifier	Limits					
13C3 HFPO-DA	91		50 - 200					

Lab Sample ID: LCSD 280-402648/3-A

Matrix: Water

Analysis Batch: 402806

Client Sample ID: Lab Control Sample Dup

Prep Type: Total/NA

Prep Batch: 402648

Analyte		Spike Added	LCSD Result	LCSD Qualifier	Unit	D	%Rec	% Rec. Limits	RPD	Limit
HFPO-DA		0.200	0.190		ug/L		95	70 - 130	1	20
Surrogate	LCSD %Recovery	LCSD Qualifier	Limits							
13C3 HFPO-DA	90		50 - 200							

Lab Sample ID: LLCS 280-402648/4-A

Matrix: Water

Analysis Batch: 402806

Client Sample ID: Lab Control Sample

Prep Type: Total/NA

Prep Batch: 402648

Analyte		Spike Added	LLCS Result	LLCS Qualifier	Unit	D	%Rec	% Rec. Limits
HFPO-DA		0.0200	0.0173		ug/L		87	70 - 130
Surrogate	LLCS %Recovery	LLCS Qualifier	Limits					
13C3 HFPO-DA	94		50 - 200					

TestAmerica Denver

# QC Sample Results

Client: Chemours Company FC, LLC The  
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 280-105698-1

## Method: 8321A - HFPO-DA (Continued)

Lab Sample ID: 280-105698-1 MS

Matrix: Water

Analysis Batch: 402806

Client Sample ID: FAY-D-6394CHKFT-W1-1-012218

Prep Type: Total/NA

Prep Batch: 402648

Analyte	Sample Result	Sample Qualifier	Spike Added	MS Result	MS Qualifier	Unit	D	%Rec	Limits
HFPO-DA	0.033		0.198	0.220		ug/L		94	70 - 130

Surrogate	MS %Recovery	MS Qualifier	Limits
13C3 HFPO-DA	86		50 - 200

Lab Sample ID: 280-105698-1 DU

Matrix: Water

Analysis Batch: 402806

Client Sample ID: FAY-D-6394CHKFT-W1-1-012218

Prep Type: Total/NA

Prep Batch: 402648

Analyte	Sample Result	Sample Qualifier	DU Result	DU Qualifier	Unit	D	RPD	Limit
HFPO-DA	0.033		0.0331		ug/L		0.9	20

Surrogate	DU %Recovery	DU Qualifier	Limits
13C3 HFPO-DA	85		50 - 200

TestAmerica Denver

# QC Association Summary

Client: Chemours Company FC, LLC The  
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 280-105698-1

## LCMS

### Analysis Batch: 390728

Lab Sample ID	Client Sample ID	Prep Type	Matrix	Method	Prep Batch
DLCK 280-390728/12	Lab Control Sample	Total/NA	Water	8321A	

### Prep Batch: 402648

Lab Sample ID	Client Sample ID	Prep Type	Matrix	Method	Prep Batch
280-105698-1	FAY-D-6394CHKFT-W1-1-012218	Total/NA	Water	3535	
280-105698-2	FAY-D-6394CHKFT-W1-1-012218-D	Total/NA	Water	3535	
280-105698-3	FAY-D-6246CHKFT-W1-1-012218	Total/NA	Water	3535	
280-105698-4	FAY-D-318BOONE-W1-1-012218	Total/NA	Water	3535	
280-105698-5	FAY-D-41BOONE-W1-1-012218	Total/NA	Water	3535	
280-105698-6	FAY-D-FB-012218	Total/NA	Water	3535	
280-105698-7	FAY-D-7145BUTLE-W1-1-012218	Total/NA	Water	3535	
280-105698-8	FAY-D-1515SCLLY-W1-1-012218	Total/NA	Water	3535	
280-105698-9	FAY-D-7396SALIE-W1-1-012218	Total/NA	Water	3535	
280-105698-10	FAY-D-7012NC87H-W1-1-012218	Total/NA	Water	3535	
MB 280-402648/1-A	Method Blank	Total/NA	Water	3535	
LCS 280-402648/2-A	Lab Control Sample	Total/NA	Water	3535	
LCSD 280-402648/3-A	Lab Control Sample Dup	Total/NA	Water	3535	
LLCS 280-402648/4-A	Lab Control Sample	Total/NA	Water	3535	
280-105698-1 MS	FAY-D-6394CHKFT-W1-1-012218	Total/NA	Water	3535	
280-105698-1 DU	FAY-D-6394CHKFT-W1-1-012218	Total/NA	Water	3535	

### Analysis Batch: 402806

Lab Sample ID	Client Sample ID	Prep Type	Matrix	Method	Prep Batch
280-105698-1	FAY-D-6394CHKFT-W1-1-012218	Total/NA	Water	8321A	402648
280-105698-2	FAY-D-6394CHKFT-W1-1-012218-D	Total/NA	Water	8321A	402648
280-105698-3	FAY-D-6246CHKFT-W1-1-012218	Total/NA	Water	8321A	402648
280-105698-4	FAY-D-318BOONE-W1-1-012218	Total/NA	Water	8321A	402648
280-105698-5	FAY-D-41BOONE-W1-1-012218	Total/NA	Water	8321A	402648
280-105698-6	FAY-D-FB-012218	Total/NA	Water	8321A	402648
280-105698-7	FAY-D-7145BUTLE-W1-1-012218	Total/NA	Water	8321A	402648
280-105698-8	FAY-D-1515SCLLY-W1-1-012218	Total/NA	Water	8321A	402648
280-105698-9	FAY-D-7396SALIE-W1-1-012218	Total/NA	Water	8321A	402648
280-105698-10	FAY-D-7012NC87H-W1-1-012218	Total/NA	Water	8321A	402648
MB 280-402648/1-A	Method Blank	Total/NA	Water	8321A	402648
LCS 280-402648/2-A	Lab Control Sample	Total/NA	Water	8321A	402648
LCSD 280-402648/3-A	Lab Control Sample Dup	Total/NA	Water	8321A	402648
LLCS 280-402648/4-A	Lab Control Sample	Total/NA	Water	8321A	402648
280-105698-1 MS	FAY-D-6394CHKFT-W1-1-012218	Total/NA	Water	8321A	402648
280-105698-1 DU	FAY-D-6394CHKFT-W1-1-012218	Total/NA	Water	8321A	402648

# Lab Chronicle

Client: Chemours Company FC, LLC The  
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 280-105698-1

Client Sample ID: FAY-D-6394CHKFT-W1-1-012218

Lab Sample ID: 280-105698-1

Date Collected: 01/22/18 16:19

Matrix: Water

Date Received: 01/23/18 09:30

Prep Type	Batch Type	Batch Method	Run	Dil Factor	Initial Amount	Final Amount	Batch Number	Prepared or Analyzed	Analyst	Lab
Total/NA	Prep	3535			263.8 mL	5 mL	402648	01/24/18 15:20	SKM	TAL DEN
Total/NA	Analysis	8321A		1			402806	01/25/18 10:54	AGCM	TAL DEN

Client Sample ID: FAY-D-6394CHKFT-W1-1-012218-D

Lab Sample ID: 280-105698-2

Date Collected: 01/22/18 16:19

Matrix: Water

Date Received: 01/23/18 09:30

Prep Type	Batch Type	Batch Method	Run	Dil Factor	Initial Amount	Final Amount	Batch Number	Prepared or Analyzed	Analyst	Lab
Total/NA	Prep	3535			250.7 mL	5 mL	402648	01/24/18 15:20	SKM	TAL DEN
Total/NA	Analysis	8321A		1			402806	01/25/18 11:07	AGCM	TAL DEN

Client Sample ID: FAY-D-6246CHKFT-W1-1-012218

Lab Sample ID: 280-105698-3

Date Collected: 01/22/18 14:56

Matrix: Water

Date Received: 01/23/18 09:30

Prep Type	Batch Type	Batch Method	Run	Dil Factor	Initial Amount	Final Amount	Batch Number	Prepared or Analyzed	Analyst	Lab
Total/NA	Prep	3535			254.1 mL	5 mL	402648	01/24/18 15:20	SKM	TAL DEN
Total/NA	Analysis	8321A		1			402806	01/25/18 11:10	AGCM	TAL DEN

Client Sample ID: FAY-D-318BOONE-W1-1-012218

Lab Sample ID: 280-105698-4

Date Collected: 01/22/18 14:15

Matrix: Water

Date Received: 01/23/18 09:30

Prep Type	Batch Type	Batch Method	Run	Dil Factor	Initial Amount	Final Amount	Batch Number	Prepared or Analyzed	Analyst	Lab
Total/NA	Prep	3535			245.1 mL	5 mL	402648	01/24/18 15:20	SKM	TAL DEN
Total/NA	Analysis	8321A		1			402806	01/25/18 11:14	AGCM	TAL DEN

Client Sample ID: FAY-D-41BOONE-W1-1-012218

Lab Sample ID: 280-105698-5

Date Collected: 01/22/18 13:49

Matrix: Water

Date Received: 01/23/18 09:30

Prep Type	Batch Type	Batch Method	Run	Dil Factor	Initial Amount	Final Amount	Batch Number	Prepared or Analyzed	Analyst	Lab
Total/NA	Prep	3535			252.7 mL	5 mL	402648	01/24/18 15:20	SKM	TAL DEN
Total/NA	Analysis	8321A		1			402806	01/25/18 11:17	AGCM	TAL DEN

Client Sample ID: FAY-D-FB-012218

Lab Sample ID: 280-105698-6

Date Collected: 01/22/18 07:30

Matrix: Water

Date Received: 01/23/18 09:30

Prep Type	Batch Type	Batch Method	Run	Dil Factor	Initial Amount	Final Amount	Batch Number	Prepared or Analyzed	Analyst	Lab
Total/NA	Prep	3535			260.4 mL	5 mL	402648	01/24/18 15:20	SKM	TAL DEN
Total/NA	Analysis	8321A		1			402806	01/25/18 11:20	AGCM	TAL DEN

TestAmerica Denver

# Lab Chronicle

Client: Chemours Company FC, LLC The  
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 280-105698-1

Client Sample ID: FAY-D-7145BUTLE-W1-1-012218

Lab Sample ID: 280-105698-7

Date Collected: 01/22/18 14:03

Matrix: Water

Date Received: 01/23/18 09:30

Prep Type	Batch Type	Batch Method	Run	Dil Factor	Initial Amount	Final Amount	Batch Number	Prepared or Analyzed	Analyst	Lab
Total/NA	Prep	3535			257.1 mL	5 mL	402648	01/24/18 15:20	SKM	TAL DEN
Total/NA	Analysis	8321A		1			402806	01/25/18 11:23	AGCM	TAL DEN

Client Sample ID: FAY-D-1515SCLLY-W1-1-012218

Lab Sample ID: 280-105698-8

Date Collected: 01/22/18 14:23

Matrix: Water

Date Received: 01/23/18 09:30

Prep Type	Batch Type	Batch Method	Run	Dil Factor	Initial Amount	Final Amount	Batch Number	Prepared or Analyzed	Analyst	Lab
Total/NA	Prep	3535			254.1 mL	5 mL	402648	01/24/18 15:20	SKM	TAL DEN
Total/NA	Analysis	8321A		1			402806	01/25/18 11:27	AGCM	TAL DEN

Client Sample ID: FAY-D-7396SALIE-W1-1-012218

Lab Sample ID: 280-105698-9

Date Collected: 01/22/18 15:39

Matrix: Water

Date Received: 01/23/18 09:30

Prep Type	Batch Type	Batch Method	Run	Dil Factor	Initial Amount	Final Amount	Batch Number	Prepared or Analyzed	Analyst	Lab
Total/NA	Prep	3535			258.6 mL	5 mL	402648	01/24/18 15:20	SKM	TAL DEN
Total/NA	Analysis	8321A		1			402806	01/25/18 11:33	AGCM	TAL DEN

Client Sample ID: FAY-D-7012NC87H-W1-1-012218

Lab Sample ID: 280-105698-10

Date Collected: 01/22/18 16:21

Matrix: Water

Date Received: 01/23/18 09:30

Prep Type	Batch Type	Batch Method	Run	Dil Factor	Initial Amount	Final Amount	Batch Number	Prepared or Analyzed	Analyst	Lab
Total/NA	Prep	3535			247.7 mL	5 mL	402648	01/24/18 15:20	SKM	TAL DEN
Total/NA	Analysis	8321A		1			402806	01/25/18 11:36	AGCM	TAL DEN

Client Sample ID: Method Blank

Lab Sample ID: MB 280-402648/1-A

Date Collected: N/A

Matrix: Water

Date Received: N/A

Prep Type	Batch Type	Batch Method	Run	Dil Factor	Initial Amount	Final Amount	Batch Number	Prepared or Analyzed	Analyst	Lab
Total/NA	Prep	3535			250 mL	5 mL	402648	01/24/18 15:20	SKM	TAL DEN
Total/NA	Analysis	8321A		1			402806	01/25/18 10:41	AGCM	TAL DEN

Client Sample ID: Lab Control Sample

Lab Sample ID: DLCK 280-390728/12

Date Collected: N/A

Matrix: Water

Date Received: N/A

Prep Type	Batch Type	Batch Method	Run	Dil Factor	Initial Amount	Final Amount	Batch Number	Prepared or Analyzed	Analyst	Lab
Total/NA	Analysis	8321A		1			390728	10/10/17 10:04	AGCM	TAL DEN

TestAmerica Denver

# Lab Chronicle

Client: Chemours Company FC, LLC The  
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 280-105698-1

Client Sample ID: Lab Control Sample

Lab Sample ID: LCS 280-402648/2-A

Date Collected: N/A

Matrix: Water

Date Received: N/A

Prep Type	Batch Type	Batch Method	Run	Dil Factor	Initial Amount	Final Amount	Batch Number	Prepared or Analyzed	Analyst	Lab
Total/NA	Prep	3535			250 mL	5 mL	402648	01/24/18 15:20	SKM	TAL DEN
Total/NA	Analysis	8321A		1			402806	01/25/18 10:44	AGCM	TAL DEN

Client Sample ID: Lab Control Sample Dup

Lab Sample ID: LCSD 280-402648/3-A

Date Collected: N/A

Matrix: Water

Date Received: N/A

Prep Type	Batch Type	Batch Method	Run	Dil Factor	Initial Amount	Final Amount	Batch Number	Prepared or Analyzed	Analyst	Lab
Total/NA	Prep	3535			250 mL	5 mL	402648	01/24/18 15:20	SKM	TAL DEN
Total/NA	Analysis	8321A		1			402806	01/25/18 10:48	AGCM	TAL DEN

Client Sample ID: Lab Control Sample

Lab Sample ID: LLCS 280-402648/4-A

Date Collected: N/A

Matrix: Water

Date Received: N/A

Prep Type	Batch Type	Batch Method	Run	Dil Factor	Initial Amount	Final Amount	Batch Number	Prepared or Analyzed	Analyst	Lab
Total/NA	Prep	3535			250 mL	5 mL	402648	01/24/18 15:20	SKM	TAL DEN
Total/NA	Analysis	8321A		1			402806	01/25/18 10:51	AGCM	TAL DEN

Client Sample ID: FAY-D-6394CHKFT-W1-1-012218

Lab Sample ID: 280-105698-1 MS

Date Collected: 01/22/18 16:19

Matrix: Water

Date Received: 01/23/18 09:30

Prep Type	Batch Type	Batch Method	Run	Dil Factor	Initial Amount	Final Amount	Batch Number	Prepared or Analyzed	Analyst	Lab
Total/NA	Prep	3535			252.1 mL	5 mL	402648	01/24/18 15:20	SKM	TAL DEN
Total/NA	Analysis	8321A		1			402806	01/25/18 11:01	AGCM	TAL DEN

## Laboratory References:

TAL DEN = TestAmerica Denver, 4955 Yarrow Street, Arvada, CO 80002, TEL (303)736-0100



## Accreditation/Certification Summary

Client: Chemours Company FC, LLC The  
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 280-105698-1

### Laboratory: TestAmerica Denver

Unless otherwise noted, all analytes for this laboratory were covered under each accreditation/certification below.

Authority	Program	EPA Region	Identification Number	Expiration Date
North Carolina (WW/SW)	State Program	4	358	12-31-18

The following analytes are included in this report, but accreditation/certification is not offered by the governing authority:

Analysis Method	Prep Method	Matrix	Analyte
8321A	3535	Water	HFPO-DA

## Method Summary

Client: Chemours Company FC, LLC The  
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 280-105698-1

---

Method	Method Description	Protocol	Laboratory
8321A	HFPO-DA	SW846	TAL DEN

---

**Protocol References:**

SW846 = "Test Methods For Evaluating Solid Waste, Physical/Chemical Methods", Third Edition, November 1986 And Its Updates.

**Laboratory References:**

TAL DEN = TestAmerica Denver, 4955 Yarrow Street, Arvada, CO 80002, TEL (303)736-0100

## Sample Summary

Client: Chemours Company FC, LLC The  
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 280-105698-1

Lab Sample ID	Client Sample ID	Matrix	Collected	Received
280-105698-1	FAY-D-6394CHKFT-W1-1-012218	Water	01/22/18 16:19	01/23/18 09:30
280-105698-2	FAY-D-6394CHKFT-W1-1-012218-D	Water	01/22/18 16:19	01/23/18 09:30
280-105698-3	FAY-D-6246CHKFT-W1-1-012218	Water	01/22/18 14:56	01/23/18 09:30
280-105698-4	FAY-D-318BOONE-W1-1-012218	Water	01/22/18 14:15	01/23/18 09:30
280-105698-5	FAY-D-41BOONE-W1-1-012218	Water	01/22/18 13:49	01/23/18 09:30
280-105698-6	FAY-D-FB-012218	Water	01/22/18 07:30	01/23/18 09:30
280-105698-7	FAY-D-7145BUTLE-W1-1-012218	Water	01/22/18 14:03	01/23/18 09:30
280-105698-8	FAY-D-1515SCLLY-W1-1-012218	Water	01/22/18 14:23	01/23/18 09:30
280-105698-9	FAY-D-7396SALIE-W1-1-012218	Water	01/22/18 15:39	01/23/18 09:30
280-105698-10	FAY-D-7012NC87H-W1-1-012218	Water	01/22/18 16:21	01/23/18 09:30

TestAmerica Denver

## LCMS MANUAL INTEGRATION SUMMARY

Lab Name: TestAmerica Denver Job No.: 280-105698-1

SDG No.: \_\_\_\_\_

Instrument ID: LC\_LCMS7 Analysis Batch Number: 390728Lab Sample ID: STD001 280-390728/3 IC Client Sample ID: \_\_\_\_\_Date Analyzed: 10/10/17 09:35 Lab File ID: hfpo717J10026.d GC Column: Synergi Hydro ID: \_\_\_\_\_

COMPOUND NAME	RETENTION TIME	MANUAL INTEGRATION		
		REASON	ANALYST	DATE
HFPO-DA	0.89	Baseline	meyera	10/10/17 11:50

## LCMS MANUAL INTEGRATION SUMMARY

Lab Name: TestAmerica Denver Job No.: 280-105698-1

SDG No.: \_\_\_\_\_

Instrument ID: LC\_LCMS7 Analysis Batch Number: 402806Lab Sample ID: LLCS 280-402648/4-A Client Sample ID: \_\_\_\_\_Date Analyzed: 01/25/18 10:51 Lab File ID: hfpo718A25010.d GC Column: Synergi Hydro ID: \_\_\_\_\_

COMPOUND NAME	RETENTION TIME	MANUAL INTEGRATION		
		REASON	ANALYST	DATE
HFPO-DA	0.96	Baseline	meyera	01/25/18 15:27

Lab Sample ID: 280-105698-1 Client Sample ID: FAY-D-6394CHKFT-W1-1-012218Date Analyzed: 01/25/18 10:54 Lab File ID: hfpo718A25011.d GC Column: Synergi Hydro ID: \_\_\_\_\_

COMPOUND NAME	RETENTION TIME	MANUAL INTEGRATION		
		REASON	ANALYST	DATE
HFPO-DA	0.96	Baseline	meyera	01/25/18 15:27

Lab Sample ID: 280-105698-2 Client Sample ID: FAY-D-6394CHKFT-W1-1-012218-DDate Analyzed: 01/25/18 11:07 Lab File ID: hfpo718A25015.d GC Column: Synergi Hydro ID: \_\_\_\_\_

COMPOUND NAME	RETENTION TIME	MANUAL INTEGRATION		
		REASON	ANALYST	DATE
HFPO-DA	1.02	Baseline	meyera	01/25/18 15:27

# REAGENT TRACEABILITY SUMMARY

Lab Name: TestAmerica Denver

Job No.: 280-105698-1

SDG No.: \_\_\_\_\_

Reagent ID	Exp Date	Prep Date	Dilutant Used	Reagent Final Volume	Parent Reagent		Analyte	Concentration
					Reagent ID	Volume Added		
<b>HFPO I.S._00007</b>	12/12/18	12/12/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA_00007	1 mL	13C3 HFPO-DA	0.5 ug/mL
							13C3 HFPO-DA (IS)	0.5 ug/mL
.13C3 HFPO-DA_00007	08/17/20	Wellington Laboratories, Lot M3HFPOADA0817			(Purchased Reagent)		13C3 HFPO-DA	50 ug/mL
							13C3 HFPO-DA (IS)	50 ug/mL
<b>HFPO Spike_00004</b>	10/30/18	10/30/17	LCMS Grade MeOH, Lot LCMS MeOH 00110	100 mL	HFPO-DA_00004	1 mL	HFPO-DA	0.5 ug/mL
.HFPO-DA_00004	07/13/20	Wellington Laboratories, Lot HFPODA0717			(Purchased Reagent)		HFPO-DA	50 ug/mL
<b>HFPO_CAL-0_00031</b>	10/24/17	10/10/17	BFC Dill_Solvent, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS MeOH 00110	100 mL	13C3 HFPO-DA_00004	1 mL	13C3 HFPO-DA	0.5 ug/mL
..13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616			(Purchased Reagent)		13C3 HFPO-DA	50 ug/mL
<b>HFPO_CAL-1_00030</b>	09/28/17	09/14/17	80:20 Methanol : H2O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L
							13C3 HFPO-DA (IS)	10 ug/L
					HFPO Spike_00003	0.5 uL	HFPO-DA	0.25 ug/L
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	13C3 HFPO-DA	0.5 ug/mL
							13C3 HFPO-DA (IS)	0.5 ug/mL
..13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616			(Purchased Reagent)		13C3 HFPO-DA	50 ug/mL
							13C3 HFPO-DA (IS)	50 ug/mL
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS MeOH 00110	100 mL	HFPO-DA_00003	1 mL	HFPO-DA	0.5 ug/mL
..HFPO-DA_00003	12/16/18	Wellington Laboratories, Lot HFPODA0213			(Purchased Reagent)		HFPO-DA	50 ug/mL
<b>HFPO_CAL-1_00031</b>	10/24/17	10/10/17	80:20 Methanol : H2O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L
							13C3 HFPO-DA (IS)	10 ug/L
					HFPO Spike_00003	0.5 uL	HFPO-DA	0.25 ug/L
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	13C3 HFPO-DA	0.5 ug/mL
							13C3 HFPO-DA (IS)	0.5 ug/mL
..13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616			(Purchased Reagent)		13C3 HFPO-DA	50 ug/mL
							13C3 HFPO-DA (IS)	50 ug/mL
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS MeOH 00110	100 mL	HFPO-DA_00003	1 mL	HFPO-DA	0.5 ug/mL
..HFPO-DA_00003	12/16/18	Wellington Laboratories, Lot HFPODA0213			(Purchased Reagent)		HFPO-DA	50 ug/mL
<b>HFPO_CAL-2_00031</b>	09/28/17	09/14/17	80:20 Methanol : H2O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L
							13C3 HFPO-DA (IS)	10 ug/L
					HFPO Spike_00003	1 uL	HFPO-DA	0.5 ug/L
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	13C3 HFPO-DA	0.5 ug/mL
							13C3 HFPO-DA (IS)	0.5 ug/mL

# REAGENT TRACEABILITY SUMMARY

Lab Name: TestAmerica Denver

Job No.: 280-105698-1

SDG No.: \_\_\_\_\_

Reagent ID	Exp Date	Prep Date	Dilutant Used	Reagent Final Volume	Parent Reagent		Analyte	Concentration
					Reagent ID	Volume Added		
..13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616			(Purchased Reagent)		13C3 HFPO-DA	50 ug/mL
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS MeOH_00110	100 mL	HFPO-DA_00003	1 mL	13C3 HFPO-DA (IS) HFPO-DA	50 ug/mL 0.5 ug/mL
..HFPO-DA_00003	12/16/18	Wellington Laboratories, Lot HFPODA0213			(Purchased Reagent)		HFPO-DA	50 ug/mL
HFPO_CAL-2_00032	10/24/17	10/10/17	80:20 Methanol : H2O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L
							13C3 HFPO-DA (IS)	10 ug/L
					HFPO Spike_00003	1 uL	HFPO-DA	0.5 ug/L
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	13C3 HFPO-DA	0.5 ug/mL
..13C3 HFPO-DA_00004	08/28/18		Wellington Laboratories, Lot M3HFPOADA0616		(Purchased Reagent)		13C3 HFPO-DA (IS)	0.5 ug/mL
							13C3 HFPO-DA	50 ug/mL
							13C3 HFPO-DA (IS)	50 ug/mL
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS MeOH_00110	100 mL	HFPO-DA_00003	1 mL	HFPO-DA	0.5 ug/mL
..HFPO-DA_00003	12/16/18	Wellington Laboratories, Lot HFPODA0213			(Purchased Reagent)		HFPO-DA	50 ug/mL
HFPO_CAL-3_00030	09/28/17	09/14/17	80:20 Methanol : H2O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L
							13C3 HFPO-DA (IS)	10 ug/L
					HFPO Spike_00003	2 uL	HFPO-DA	1 ug/L
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	13C3 HFPO-DA	0.5 ug/mL
..13C3 HFPO-DA_00004	08/28/18		Wellington Laboratories, Lot M3HFPOADA0616		(Purchased Reagent)		13C3 HFPO-DA (IS)	0.5 ug/mL
							13C3 HFPO-DA	50 ug/mL
							13C3 HFPO-DA (IS)	50 ug/mL
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS MeOH_00110	100 mL	HFPO-DA_00003	1 mL	HFPO-DA	0.5 ug/mL
..HFPO-DA_00003	12/16/18	Wellington Laboratories, Lot HFPODA0213			(Purchased Reagent)		HFPO-DA	50 ug/mL
HFPO_CAL-3_00031	10/24/17	10/10/17	80:20 Methanol : H2O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L
							13C3 HFPO-DA (IS)	10 ug/L
					HFPO Spike_00003	2 uL	HFPO-DA	1 ug/L
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	13C3 HFPO-DA	0.5 ug/mL
..13C3 HFPO-DA_00004	08/28/18		Wellington Laboratories, Lot M3HFPOADA0616		(Purchased Reagent)		13C3 HFPO-DA (IS)	0.5 ug/mL
							13C3 HFPO-DA	50 ug/mL
							13C3 HFPO-DA (IS)	50 ug/mL
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS MeOH_00110	100 mL	HFPO-DA_00003	1 mL	HFPO-DA	0.5 ug/mL
..HFPO-DA_00003	12/16/18	Wellington Laboratories, Lot HFPODA0213			(Purchased Reagent)		HFPO-DA	50 ug/mL
HFPO_CAL-4_00030	09/28/17	09/14/17	80:20 Methanol : H2O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L
							13C3 HFPO-DA (IS)	10 ug/L

# REAGENT TRACEABILITY SUMMARY

Lab Name: TestAmerica Denver

Job No.: 280-105698-1

SDG No.: \_\_\_\_\_

Reagent ID	Exp Date	Prep Date	Dilutant Used	Reagent Final Volume	Parent Reagent		Analyte	Concentration
					Reagent ID	Volume Added		
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	HFPO Spike 00003	4 uL	HFPO-DA	2 ug/L
..13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616			13C3 HFPO-DA_00004	1 mL	13C3 HFPO-DA	0.5 ug/mL
					(Purchased Reagent)		13C3 HFPO-DA (IS)	0.5 ug/mL
							13C3 HFPO-DA	50 ug/mL
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS MeOH_00110	100 mL	HFPO-DA_00003	1 mL	13C3 HFPO-DA (IS)	50 ug/mL
..HFPO-DA_00003	12/16/18	Wellington Laboratories, Lot HFPODA0213			(Purchased Reagent)		HFPO-DA	0.5 ug/mL
<b>HFPO_CAL-4_00031</b>	10/24/17	10/10/17	80:20 Methanol : H2O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	HFPO-DA	50 ug/mL
							13C3 HFPO-DA	10 ug/L
					HFPO Spike 00003	4 uL	13C3 HFPO-DA (IS)	10 ug/L
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	HFPO-DA	2 ug/L
..13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616			(Purchased Reagent)		13C3 HFPO-DA	0.5 ug/mL
							13C3 HFPO-DA (IS)	0.5 ug/mL
							13C3 HFPO-DA	50 ug/mL
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS MeOH_00110	100 mL	HFPO-DA_00003	1 mL	13C3 HFPO-DA (IS)	50 ug/mL
..HFPO-DA_00003	12/16/18	Wellington Laboratories, Lot HFPODA0213			(Purchased Reagent)		HFPO-DA	0.5 ug/mL
<b>HFPO_CAL-5_00067</b>	09/28/17	09/14/17	80:20 Methanol : H2O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L
							13C3 HFPO-DA (IS)	10 ug/L
					HFPO Spike 00003	10 uL	HFPO-DA	5 ug/L
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	13C3 HFPO-DA	0.5 ug/mL
..13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616			(Purchased Reagent)		13C3 HFPO-DA (IS)	0.5 ug/mL
							13C3 HFPO-DA	50 ug/mL
							13C3 HFPO-DA (IS)	50 ug/mL
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS MeOH_00110	100 mL	HFPO-DA_00003	1 mL	HFPO-DA	0.5 ug/mL
..HFPO-DA_00003	12/16/18	Wellington Laboratories, Lot HFPODA0213			(Purchased Reagent)		HFPO-DA	50 ug/mL
<b>HFPO_CAL-5_00070</b>	10/24/17	10/10/17	80:20 Methanol : H2O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L
							13C3 HFPO-DA (IS)	10 ug/L
					HFPO Spike 00003	10 uL	HFPO-DA	5 ug/L
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	13C3 HFPO-DA	0.5 ug/mL
..13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616			(Purchased Reagent)		13C3 HFPO-DA (IS)	0.5 ug/mL
							13C3 HFPO-DA	50 ug/mL
							13C3 HFPO-DA (IS)	50 ug/mL
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS MeOH_00110	100 mL	HFPO-DA_00003	1 mL	HFPO-DA	0.5 ug/mL
..HFPO-DA_00003	12/16/18	Wellington Laboratories, Lot HFPODA0213			(Purchased Reagent)		HFPO-DA	50 ug/mL



# REAGENT TRACEABILITY SUMMARY

Lab Name: TestAmerica Denver

Job No.: 280-105698-1

SDG No.: \_\_\_\_\_

Reagent ID	Exp Date	Prep Date	Dilutant Used	Reagent Final Volume	Parent Reagent		Analyte	Concentration
					Reagent ID	Volume Added		
<b>HFPO_CAL-5_00078</b>	01/26/18	01/12/18	80:20 Methanol : H2O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L
					HFPO Spike_00004	10 uL	HFPO-DA	5 ug/L
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	13C3 HFPO-DA	0.5 ug/mL
..13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616			(Purchased Reagent)		13C3 HFPO-DA	50 ug/mL
.HFPO Spike_00004	10/30/18	10/30/17	LCMS Grade MeOH, Lot LCMS MeOH_00110	100 mL	HFPO-DA_00004	1 mL	HFPO-DA	0.5 ug/mL
..HFPO-DA_00004	07/13/20	Wellington Laboratories, Lot HFPODA0717			(Purchased Reagent)		HFPO-DA	50 ug/mL
<b>HFPO_CAL-6_00067</b>	09/28/17	09/14/17	80:20 Methanol : H2O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L
							13C3 HFPO-DA (IS)	10 ug/L
					HFPO Spike_00003	20 uL	HFPO-DA	10 ug/L
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	13C3 HFPO-DA	0.5 ug/mL
							13C3 HFPO-DA (IS)	0.5 ug/mL
..13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616			(Purchased Reagent)		13C3 HFPO-DA	50 ug/mL
							13C3 HFPO-DA (IS)	50 ug/mL
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS MeOH_00110	100 mL	HFPO-DA_00003	1 mL	HFPO-DA	0.5 ug/mL
..HFPO-DA_00003	12/16/18	Wellington Laboratories, Lot HFPODA0213			(Purchased Reagent)		HFPO-DA	50 ug/mL
<b>HFPO_CAL-6_00070</b>	10/24/17	10/10/17	80:20 Methanol : H2O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L
							13C3 HFPO-DA (IS)	10 ug/L
					HFPO Spike_00003	20 uL	HFPO-DA	10 ug/L
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	13C3 HFPO-DA	0.5 ug/mL
							13C3 HFPO-DA (IS)	0.5 ug/mL
..13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616			(Purchased Reagent)		13C3 HFPO-DA	50 ug/mL
							13C3 HFPO-DA (IS)	50 ug/mL
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS MeOH_00110	100 mL	HFPO-DA_00003	1 mL	HFPO-DA	0.5 ug/mL
..HFPO-DA_00003	12/16/18	Wellington Laboratories, Lot HFPODA0213			(Purchased Reagent)		HFPO-DA	50 ug/mL
<b>HFPO_CAL-6_00078</b>	01/26/18	01/12/18	80:20 Methanol : H2O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L
					HFPO Spike_00004	20 uL	HFPO-DA	10 ug/L
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	13C3 HFPO-DA	0.5 ug/mL
..13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616			(Purchased Reagent)		13C3 HFPO-DA	50 ug/mL
.HFPO Spike_00004	10/30/18	10/30/17	LCMS Grade MeOH, Lot LCMS MeOH_00110	100 mL	HFPO-DA_00004	1 mL	HFPO-DA	0.5 ug/mL
..HFPO-DA_00004	07/13/20	Wellington Laboratories, Lot HFPODA0717			(Purchased Reagent)		HFPO-DA	50 ug/mL
<b>HFPO_CAL-7_00030</b>	09/28/17	09/14/17	80:20 Methanol : H2O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L

# REAGENT TRACEABILITY SUMMARY

Lab Name: TestAmerica Denver

Job No.: 280-105698-1

SDG No.: \_\_\_\_\_

Reagent ID	Exp Date	Prep Date	Dilutant Used	Reagent Final Volume	Parent Reagent		Analyte	Concentration
					Reagent ID	Volume Added		
							13C3 HFPO-DA (IS)	10 ug/L
					HFPO Spike_00003	50 uL	HFPO-DA	25 ug/L
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA__00004	1 mL	13C3 HFPO-DA	0.5 ug/mL
..13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616			(Purchased Reagent)		13C3 HFPO-DA (IS)	0.5 ug/mL
							13C3 HFPO-DA	50 ug/mL
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS MeOH_00110	100 mL	HFPO-DA_00003	1 mL	13C3 HFPO-DA (IS)	50 ug/mL
							HFPO-DA	0.5 ug/mL
..HFPO-DA_00003	12/16/18	Wellington Laboratories, Lot HFPODA0213			(Purchased Reagent)		HFPO-DA	50 ug/mL
<b>HFPO_CAL-7_00031</b>	10/24/17	10/10/17	80:20 Methanol : H2O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L
							13C3 HFPO-DA (IS)	10 ug/L
					HFPO Spike_00003	50 uL	HFPO-DA	25 ug/L
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA__00004	1 mL	13C3 HFPO-DA	0.5 ug/mL
..13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616			(Purchased Reagent)		13C3 HFPO-DA (IS)	0.5 ug/mL
							13C3 HFPO-DA	50 ug/mL
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS MeOH_00110	100 mL	HFPO-DA_00003	1 mL	13C3 HFPO-DA (IS)	50 ug/mL
							HFPO-DA	0.5 ug/mL
..HFPO-DA_00003	12/16/18	Wellington Laboratories, Lot HFPODA0213			(Purchased Reagent)		HFPO-DA	50 ug/mL
<b>HFPO_CAL-8_00030</b>	09/28/17	09/14/17	80:20 Methanol : H2O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L
							13C3 HFPO-DA (IS)	10 ug/L
					HFPO Spike_00003	100 uL	HFPO-DA	50 ug/L
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA__00004	1 mL	13C3 HFPO-DA	0.5 ug/mL
..13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616			(Purchased Reagent)		13C3 HFPO-DA (IS)	0.5 ug/mL
							13C3 HFPO-DA	50 ug/mL
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS MeOH_00110	100 mL	HFPO-DA_00003	1 mL	13C3 HFPO-DA (IS)	50 ug/mL
							HFPO-DA	0.5 ug/mL
..HFPO-DA_00003	12/16/18	Wellington Laboratories, Lot HFPODA0213			(Purchased Reagent)		HFPO-DA	50 ug/mL
<b>HFPO_CAL-8_00031</b>	10/24/17	10/10/17	80:20 Methanol : H2O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L
							13C3 HFPO-DA (IS)	10 ug/L
					HFPO Spike_00003	100 uL	HFPO-DA	50 ug/L
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA__00004	1 mL	13C3 HFPO-DA	0.5 ug/mL
..13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616			(Purchased Reagent)		13C3 HFPO-DA (IS)	0.5 ug/mL
							13C3 HFPO-DA	50 ug/mL
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS MeOH_00110	100 mL	HFPO-DA_00003	1 mL	13C3 HFPO-DA (IS)	50 ug/mL
							HFPO-DA	0.5 ug/mL

REAGENT TRACEABILITY SUMMARY

Lab Name: TestAmerica Denver

Job No.: 280-105698-1


SDG No.: \_\_\_\_\_

Reagent ID	Exp Date	Prep Date	Dilutant Used	Reagent Final Volume	Parent Reagent		Analyte	Concentration
					Reagent ID	Volume Added		
..HFPO-DA_00003	12/16/18	Wellington Laboratories, Lot HFPODA0213			(Purchased Reagent)		HFPO-DA	50 ug/mL
HFPO_ICV_00031	09/28/17	09/14/17	80:20 Methanol : H2O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L
					HFPO Spike 00003	4 uL	HFPO-DA	2 ug/L
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS MeOH 00110	100 mL	13C3 HFPO-DA_00004	1 mL	13C3 HFPO-DA	0.5 ug/mL
..13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616			(Purchased Reagent)		13C3 HFPO-DA	50 ug/mL
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS MeOH 00110	100 mL	HFPO-DA_00003	1 mL	HFPO-DA	0.5 ug/mL
..HFPO-DA_00003	12/16/18	Wellington Laboratories, Lot HFPODA0213			(Purchased Reagent)		HFPO-DA	50 ug/mL
HFPO_ICV_00032	10/24/17	10/10/17	80:20 Methanol : H2O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L
					HFPO Spike_00003	4 uL	HFPO-DA	2 ug/L
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS MeOH 00110	100 mL	13C3 HFPO-DA_00004	1 mL	13C3 HFPO-DA	0.5 ug/mL
..13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616			(Purchased Reagent)		13C3 HFPO-DA	50 ug/mL
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS MeOH 00110	100 mL	HFPO-DA_00003	1 mL	HFPO-DA	0.5 ug/mL
..HFPO-DA_00003	12/16/18	Wellington Laboratories, Lot HFPODA0213			(Purchased Reagent)		HFPO-DA	50 ug/mL

Reagent

---

**13C3 HFPO-DA\_00004**

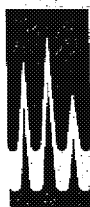

# WELLINGTON LABORATORIES

## CERTIFICATE OF ANALYSIS DOCUMENTATION

**PRODUCT CODE:**


M3HFPO-DA

**LOT NUMBER:**

M3HFPODA0616

**COMPOUND:**2,3,3,3-Tetrafluoro-2-(1,1,2,2,3,3,3-heptafluoropropoxy)-<sup>13</sup>C<sub>3</sub>-propanoic acid**STRUCTURE:****CAS #:**

Not available

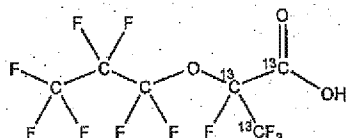**MOLECULAR FORMULA:** $^{13}\text{C}_3^{12}\text{C}_3\text{HF}_{11}\text{O}_3$ **CONCENTRATION:**

50 ± 2.5 µg/ml

**CHEMICAL PURITY:**

&gt;98%

**LAST TESTED:** (mm/dd/yyyy)

06/25/2016

**EXPIRY DATE:** (mm/dd/yyyy)

06/25/2019

**RECOMMENDED STORAGE:**

Store ampoule in a cool, dark place

**MOLECULAR WEIGHT:**

333.03

**SOLVENT(S):**

Methanol

**ISOTOPIC PURITY:**≥99% <sup>13</sup>C  
(<sup>13</sup>C<sub>3</sub>)**DOCUMENTATION/ DATA ATTACHED:**

Figure 1: LC/MS Data (TIC and Mass Spectrum)

Figure 2: LC/MS/MS Data (Selected MRM Transitions)

**ADDITIONAL INFORMATION:**

- See page 2 for further details.
- Contains ~ 1.5% of two constitutional isomers.
- Product is commercially known as GenX.

**FOR LABORATORY USE ONLY: NOT FOR HUMAN OR DRUG USE**

Certified By:

  
B.G. Chittim

Date:

06/29/2016

(mm/dd/yyyy)

Wellington Laboratories Inc., 345 Southgate Dr. Guelph ON N1G 3M5 CANADA  
519-822-2436 • Fax: 519-822-2849 • info@well-labs.com

#### **INTENDED USE:**

The products prepared by Wellington Laboratories Inc. are for laboratory use only. This certified reference material (CRM) was designed to be used as a standard for the identification and/or quantification of the specific chemical compound it contains.

#### **HAZARDS:**

This product should only be used by qualified personnel familiar with its potential hazards and trained in the handling of hazardous chemicals. Due care should be exercised to prevent unnecessary human contact or ingestion. All procedures should be carried out in a well-functioning fume hood and suitable gloves, eye protection, and clothing should be worn at all times. Waste should be disposed of according to national and regional regulations. Safety Data Sheets (SDSs) are available upon request.

#### **SYNTHESIS / CHARACTERIZATION:**

Where possible, all of our products are synthesized using single-product unambiguous routes. They are then characterized, and their structures and purities confirmed, using a combination of the most relevant techniques, such as NMR, GC/MS, LC/MS/MS, SFC/UV/MS/MS, x-ray crystallography, and melting point. Isotopic purities of mass-labelled compounds are also confirmed using HRGC/HRMS and/or LC/MS/MS.

#### **HOMOGENEITY:**

Prior to solution preparation, crystalline material is tested for homogeneity using a variety of techniques (as stated above) and its solubility in a given diluent is taken into consideration. Duplicate solutions of a new product are prepared from the same crystalline lot and, after the addition of an appropriate internal standard, they are compared by GC/MS, LC/MS/MS and/or SFC/UV/MS/MS. The relative response factors of the analyte of interest in each solution are required to be <5% RSD. New solution lots of existing products are compared to older lots in the same manner, which further confirms the homogeneity of the crystalline material as well as the stability and homogeneity of the solutions in the storage containers.

#### **UNCERTAINTY:**

The maximum combined relative standard uncertainty of our reference standard solutions is calculated using the following equation:

The combined relative standard uncertainty,  $u_c(y)$ , of a value  $y$  and the uncertainty of the independent parameters  $x_1, x_2, \dots, x_n$  on which it depends is:

$$u_c(y(x_1, x_2, \dots, x_n)) = \sqrt{\sum_{i=1}^n u(y, x_i)^2}$$

where  $x$  is expressed as a relative standard uncertainty of the individual parameter.

The individual uncertainties taken into account include those associated with weights (calibration of the balance) and volumes (calibration of the volumetric glassware). An expanded maximum combined percent relative uncertainty of  $\pm 5\%$  (calculated with a coverage factor of 2 and a level of confidence of 95%) is stated on the Certificate of Analysis for all of our products.

#### **TRACEABILITY:**

All reference standard solutions are traceable to specific crystalline lots. The microbalances used for solution preparation are regularly tested by an external ISO/IEC 17025 accredited calibration company. In addition, their calibration is verified prior to each weighing using NIST and/or NRC traceable external weights. All volumetric glassware used is of Class A tolerance and has been tested according to the appropriate ASTM procedures, which are ultimately traceable to NIST. For certain products, traceability to international interlaboratory studies has also been established.

#### **EXPIRY DATE / PERIOD OF VALIDITY:**

Ongoing stability studies of this product have demonstrated stability in its composition and concentration, until the specified expiry date, in the unopened ampoule. Monitoring for any degradation or change in concentration of the listed analyte(s) is performed on a routine basis.

#### **LIMITED WARRANTY:**

At the time of shipment, all products are warranted to be free of defects in material and workmanship and to conform to the stated technical and purity specifications.



#### **QUALITY MANAGEMENT:**

This product was produced using a Quality Management System registered to the latest versions of ISO 9001 by SAI Global, ISO/IEC 17025 by the Canadian Association for Laboratory Accreditation Inc. (CALA; A 1226), and ISO GUIDE 34 by ANSI-ASQ National Accreditation Board (ANAB; AR-1523).



\*\*For additional information or assistance concerning this or any other products from Wellington Laboratories Inc., please visit our website at [www.well-labs.com](http://www.well-labs.com) or contact us directly at [info@well-labs.com](mailto:info@well-labs.com)\*\*

**Figure 1: M3HFPO-DA; LC/MS Data (TIC and Mass Spectrum)**

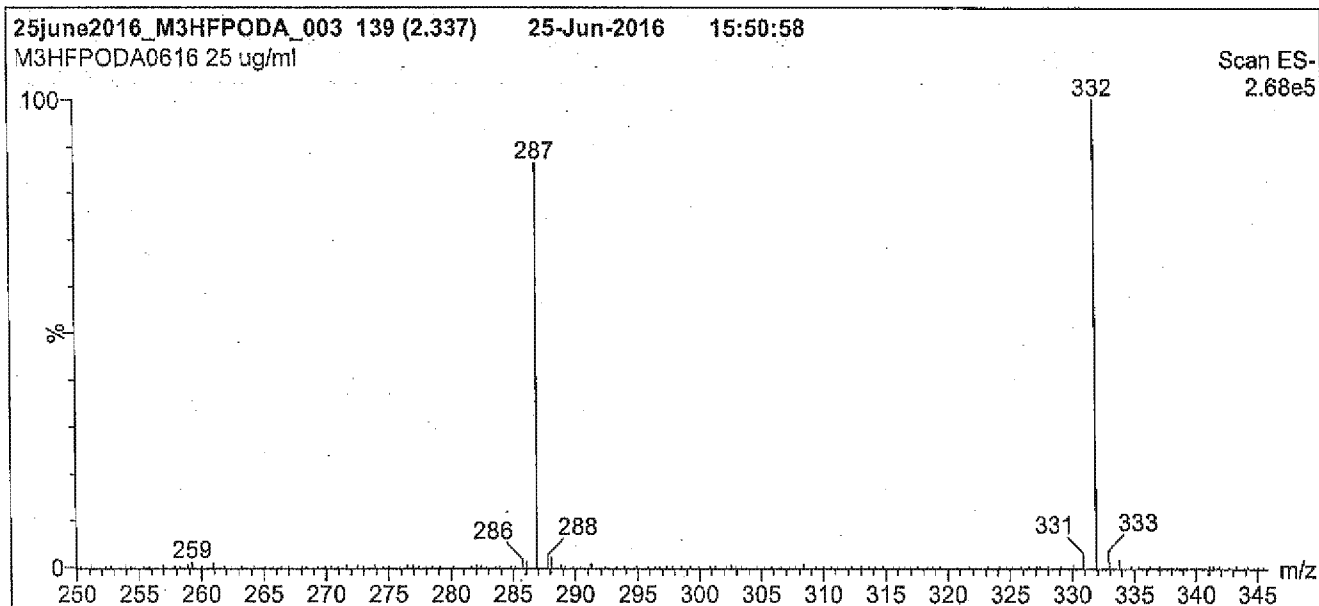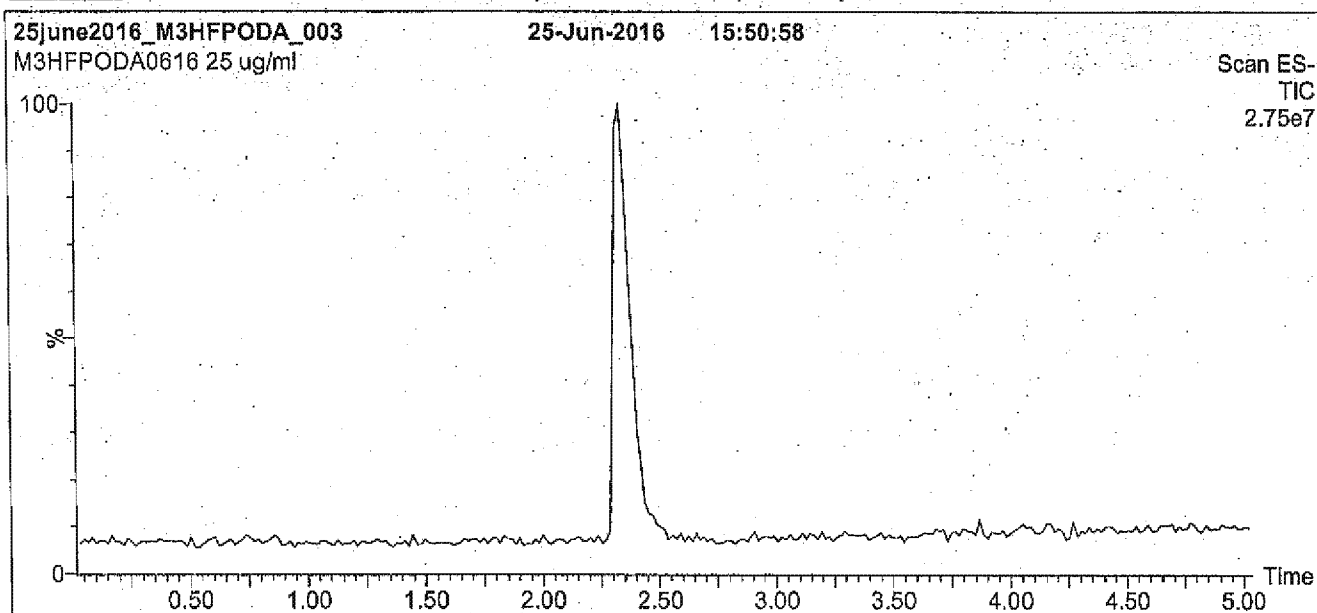

**Conditions for Figure 1:**

**LC:** Waters Acquity Ultra Performance LC  
**MS:** Micromass Quattro micro API MS

**Chromatographic Conditions**

Column: Acquity UPLC BEH Shield RP<sub>18</sub>  
1.7  $\mu$ m, 2.1 x 100 mm

Mobile phase: Gradient

Start: 50% (80:20 MeOH:ACN) / 50% H<sub>2</sub>O  
(both with 10 mM NH<sub>4</sub>OAc buffer)  
Ramp to 90% organic over 7 min and hold for 1.5 min  
before returning to initial conditions in 0.5 min.  
Time: 10 min


Flow: 300  $\mu$ l/min

**MS Parameters**

Experiment: Full Scan (250 - 850 amu)

Source: Electrospray (negative)  
Capillary Voltage (kV) = 3.00  
Cone Voltage (V) = 9.00  
Cone Gas Flow (l/hr) = 50  
Desolvation Gas Flow (l/hr) = 750

**Figure 2:** M3HFPO-DA; LC/MS/MS Data (Selected MRM Transitions)

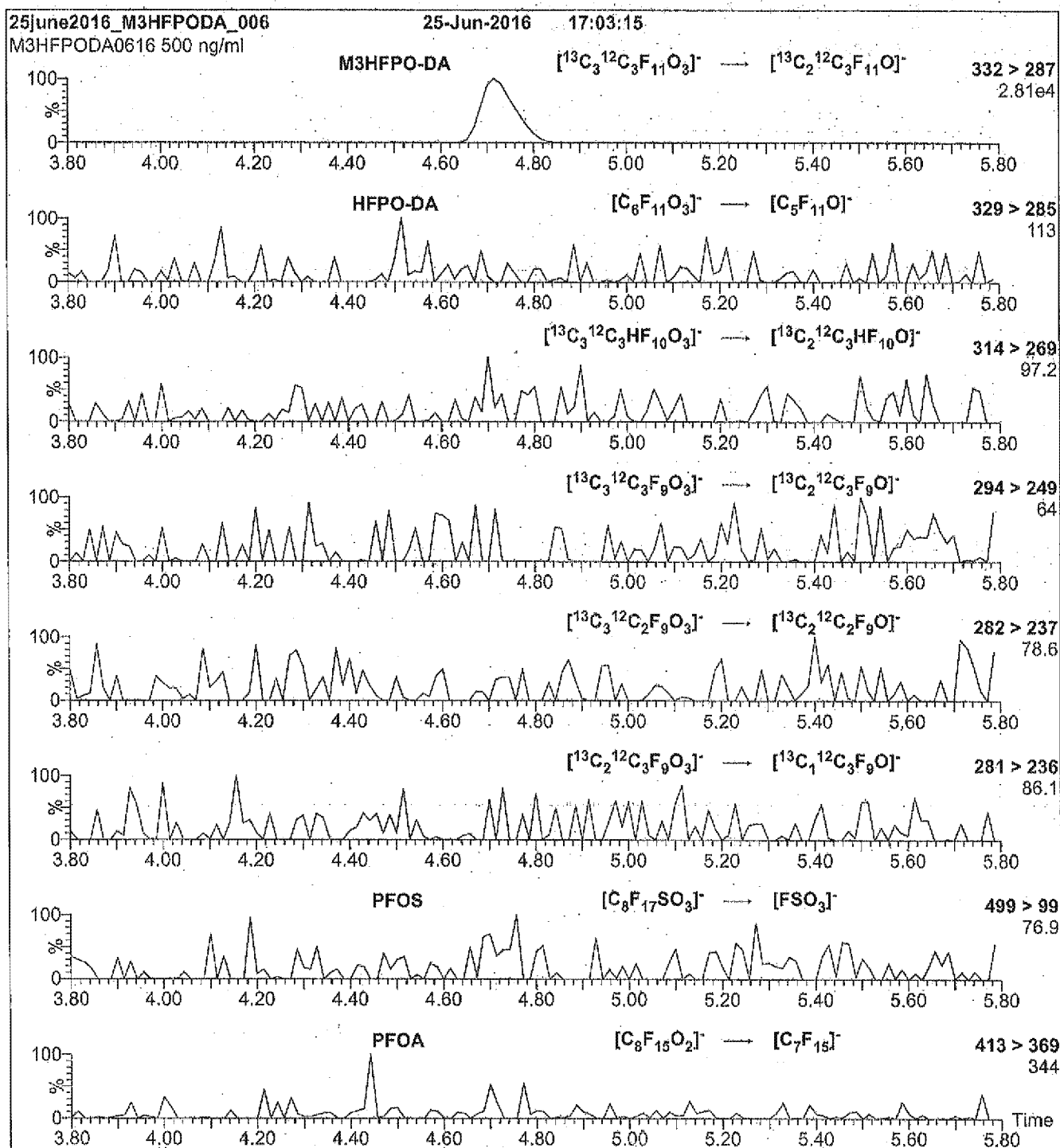

**Conditions for Figure 2:**

Injection: Direct loop Injection  
10  $\mu\text{l}$  (500 ng/ml M3HFPO-DA)

Mobile phase: Isocratic 80% (80:20 MeOH:ACN) / 20%  $\text{H}_2\text{O}$   
(both with 10 mM  $\text{NH}_4\text{OAc}$  buffer)

Flow: 300  $\mu\text{l}/\text{min}$

**MS Parameters**

Collision Gas (mbar) =  $3.46 \times 10^{-3}$   
Collision Energy (eV) = 5



Reagent

---

**13C3 HFPO-DA\_00007**



# WELLINGTON LABORATORIES


## CERTIFICATE OF ANALYSIS DOCUMENTATION

**PRODUCT CODE:**

M3HFPO-DA

**LOT NUMBER:** M3HFPODA0817**COMPOUND:**2,3,3,3-Tetrafluoro-2-(1,1,2,2,3,3,3-heptafluoropropoxy)-<sup>13</sup>C<sub>3</sub>-propanoic acid**STRUCTURE:****CAS #:**

Not available

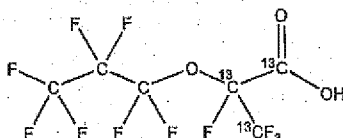**MOLECULAR FORMULA:**<sup>13</sup>C<sub>3</sub><sup>12</sup>C<sub>3</sub>HF<sub>11</sub>O<sub>3</sub>**CONCENTRATION:**

50 ± 2.5 µg/ml

**CHEMICAL PURITY:**

&gt;98%

**LAST TESTED:** (mm/dd/yyyy)

08/17/2017

**EXPIRY DATE:** (mm/dd/yyyy)

08/17/2020

**RECOMMENDED STORAGE:**

Store ampoule in a cool, dark place

**MOLECULAR WEIGHT:**

333.03

**SOLVENT(S):**

Methanol

**ISOTOPIC PURITY:**≥99% <sup>13</sup>C(<sup>13</sup>C<sub>3</sub>)**DOCUMENTATION DATA ATTACHED:**

Figure 1: LC/MS Data (TIC and Mass Spectrum)

Figure 2: LC/MS/MS Data (Selected MRM Transitions)

**ADDITIONAL INFORMATION:**

- See page 2 for further details.
- Contains ~ 1.5% of two constitutional isomers.
- Product is commercially known as GenX.

**FOR LABORATORY USE ONLY: NOT FOR HUMAN OR DRUG USE**

Certified By:

  
B.G. Chittim, General Manager

Date: 08/25/2017

(mm/dd/yyyy)

Wellington Laboratories Inc., 345 Southgate Dr. Guelph ON N1G 3M5 CANADA  
519-822-2436 • Fax: 519-822-2849 • Info@well-labs.com

### **INTENDED USE:**

The products prepared by Wellington Laboratories Inc. are for laboratory use only. This certified reference material (CRM) was designed to be used as a standard for the identification and/or quantification of the specific chemical compound it contains.

### **HAZARDS:**

This product should only be used by qualified personnel familiar with its potential hazards and trained in the handling of hazardous chemicals. Due care should be exercised to prevent unnecessary human contact or ingestion. All procedures should be carried out in a well-functioning fume hood and suitable gloves, eye protection, and clothing should be worn at all times. Waste should be disposed of according to national and regional regulations. Safety Data Sheets (SDSs) are available upon request.

### **SYNTHESIS / CHARACTERIZATION:**

Where possible, all of our products are synthesized using single-product unambiguous routes. They are then characterized, and their structures and purities confirmed, using a combination of the most relevant techniques, such as NMR, GC/MS, LC/MS/MS, SFC/UV/MS/MS, x-ray crystallography, and melting point. Isotopic purities of mass-labelled compounds are also confirmed using HRGC/HRMS and/or LC/MS/MS.

### **HOMOGENEITY:**

Prior to solution preparation, crystalline material is tested for homogeneity using a variety of techniques (as stated above) and its solubility in a given diluent is taken into consideration. Duplicate solutions of a new product are prepared from the same crystalline lot and, after the addition of an appropriate internal standard, they are compared by GC/MS, LC/MS/MS and/or SFC/UV/MS/MS. The relative response factors of the analyte of interest in each solution are required to be <5% RSD. New solution lots of existing products are compared to older lots in the same manner, which further confirms the homogeneity of the crystalline material as well as the stability and homogeneity of the solutions in the storage containers. In order to maintain the integrity of the assigned value(s), and associated uncertainty, the dilution or injection of a subsample of this product should be performed using calibrated measuring equipment.

### **UNCERTAINTY:**

The maximum combined relative standard uncertainty of our reference standard solutions is calculated using the following equation:

The combined relative standard uncertainty,  $u_c(y)$ , of a value  $y$  and the uncertainty of the independent parameters

$x_1, x_2, \dots, x_n$  on which it depends is:

$$u_c(y(x_1, x_2, \dots, x_n)) = \sqrt{\sum_{i=1}^n u(y, x_i)^2}$$

where  $x$  is expressed as a relative standard uncertainty of the individual parameter.

The individual uncertainties taken into account include those associated with weights (calibration of the balance) and volumes (calibration of the volumetric glassware). An expanded maximum combined percent relative uncertainty of  $\pm 5\%$  (calculated with a coverage factor of 2 and a level of confidence of 95%) is stated on the Certificate of Analysis for all of our products.

### **TRACEABILITY:**

All reference standard solutions are traceable to specific crystalline lots. The microbalances used for solution preparation are regularly tested by an external ISO/IEC 17025 accredited calibration company. In addition, their calibration is verified prior to each weighing using calibrated NIST and/or NRC traceable external weights. All volumetric glassware used is calibrated, of Class A tolerance, and has been tested according to the appropriate ASTM procedures, which are ultimately traceable to NIST. For certain products, traceability to international interlaboratory studies has also been established.

### **EXPIRY DATE / PERIOD OF VALIDITY:**

Ongoing stability studies of this product have demonstrated stability in its composition and concentration, until the specified expiry date, in the unopened ampoule. Monitoring for any degradation or change in concentration of the listed analyte(s) is performed on a routine basis.

### **LIMITED WARRANTY:**

At the time of shipment, all products are warranted to be free of defects in material and workmanship and to conform to the stated technical and purity specifications.



### **QUALITY MANAGEMENT:**

This product was produced using a Quality Management System registered to the latest versions of ISO 9001 by SAI Global, ISO/IEC 17025 by the Canadian Association for Laboratory Accreditation Inc. (CALA; A 1226), and ISO GUIDE 34 by ANSI-ASQ National Accreditation Board (ANAB; AR-1523).



\*\*For additional information or assistance concerning this or any other products from Wellington Laboratories Inc., please visit our website at [www.well-labs.com](http://www.well-labs.com) or contact us directly at [info@well-labs.com](mailto:info@well-labs.com)\*\*

**Figure 1: M3HFPO-DA; LC/MS Data (TIC and Mass Spectrum)**

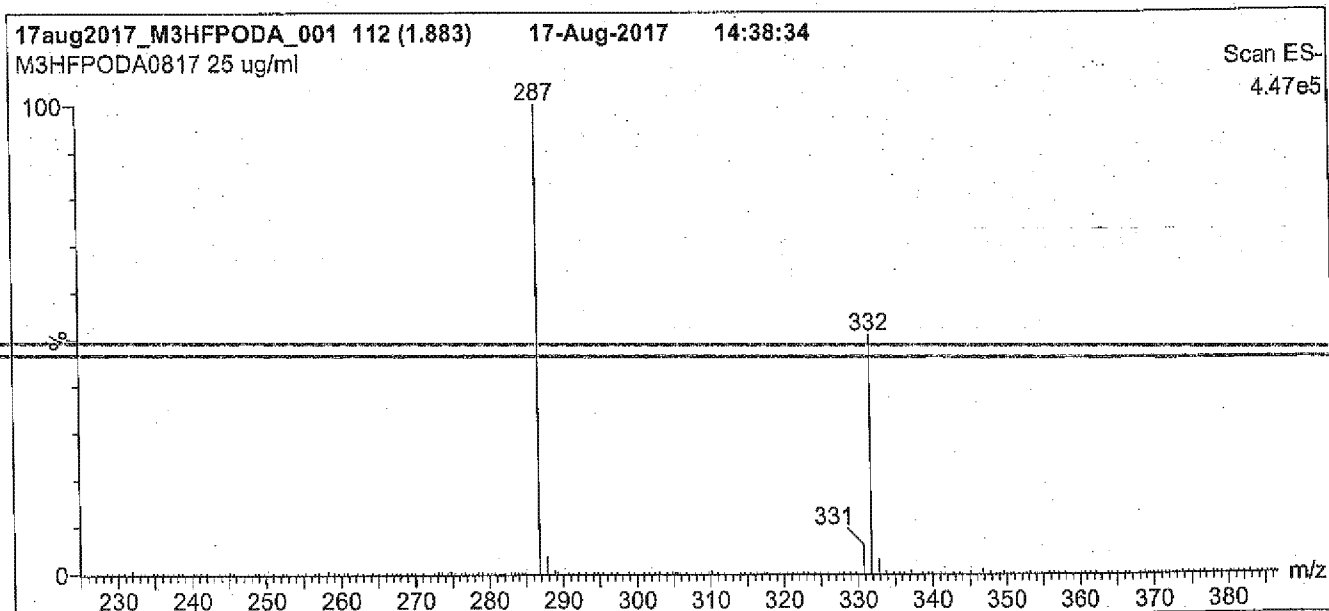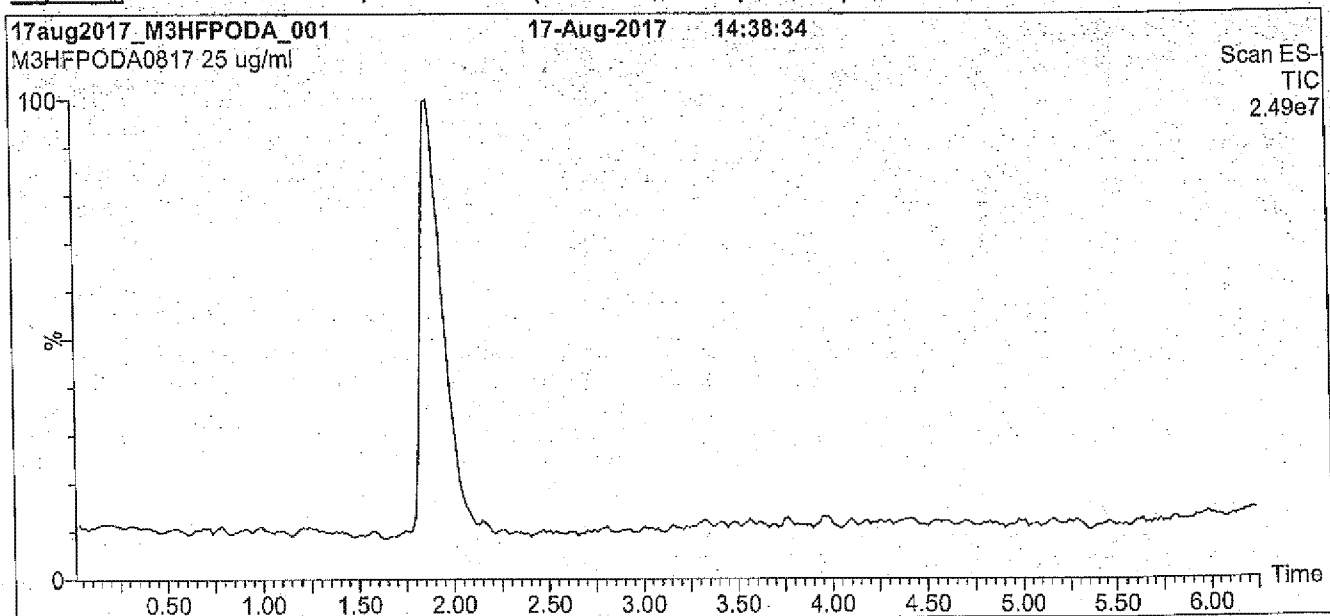

**Conditions for Figure 1:**

**LC:** Waters Acquity Ultra Performance LC  
**MS:** Micromass Quattro *micro* API MS

**Chromatographic Conditions**

**Column:** Acquity UPLC BEH Shield RP<sub>18</sub>  
1.7  $\mu$ m, 2.1 x 100 mm

**Mobile phase:** Gradient  
Start: 55% MeOH / 45% H<sub>2</sub>O with 10 mM NH<sub>4</sub>OAc buffer  
Ramp to 90% organic over 7.5 min and hold for 1.5 min  
before returning to initial conditions in 0.5 min.

**Time:** 10 min


**Flow:** 300  $\mu$ l/min

**MS Parameters**

**Experiment:** Full Scan (225 - 850 amu)

**Source:** Electrospray (negative)  
**Capillary Voltage (kV)** = 3.00  
**Cone Voltage (V)** = 10.00  
**Cone Gas Flow (l/hr)** = 100  
**Desolvation Gas Flow (l/hr)** = 750

**Figure 2: M3HFPO-DA; LC/MS/MS Data (Selected MRM Transitions)**

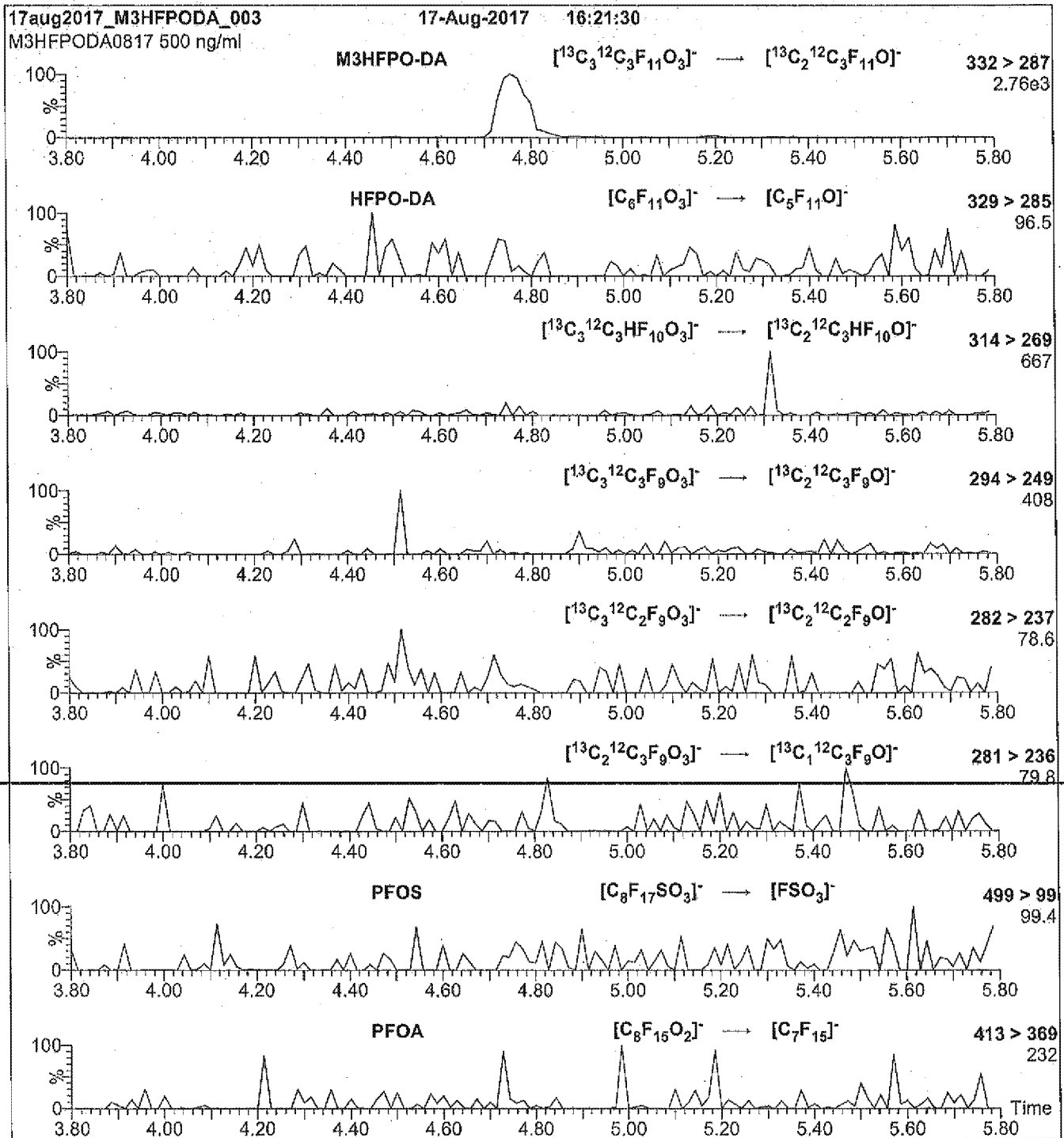

**Conditions for Figure 2:**

Injection: Direct loop Injection  
10  $\mu\text{l}$  (500 ng/ml M3HFPO-DA)

Mobile phase: Isocratic 80% MeOH / 20%  $\text{H}_2\text{O}$  with 10 mM  $\text{NH}_4\text{OAc}$  buffer

Flow: 300  $\mu\text{l}/\text{min}$


**MS Parameters**

Collision Gas (mbar) =  $3.24\text{e-}3$   
Collision Energy (eV) = 5

Reagent

---

**HFPO I.S.\_00004**

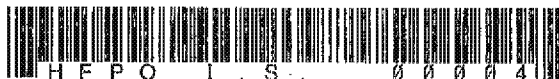

**Reagent ID:** HFPO I.S.\_00004

Description: Internal Standard for HFPO 0.5ug/ml  
 No. of Bottles: 1  
 Storage Location: North Analytical  
 Reagent Volume: 100.000 mL  
 Creation Date: 08/28/2017  
 Open Date:  
 Container(s): 4700620  
 Comment:

Expiration Date: 08/28/2018  
 Laboratory: TestAmerica Denver  
 Prepared By: Meyer, Andrew GC  
 Solvent: LCMS Grade MeOH  
 Solvent Lot: LCMS\_MeOH\_00110

## Reagent Analyte Information

Analyte	Source ID	Source Exp. Date	Source Conc.	Source Conc. Units	Final Conc.	Final Conc. Units
13C3 HFPO-DA	13C3 HFPO-DA_00004	08/28/2018	50.00000	ug/mL	0.50000	ug/mL
13C3 HFPO-DA (IS)	13C3 HFPO-DA_00004	08/28/2018	50.00000	ug/mL	0.50000	ug/mL

## Source Reagents

Reagent	Description	Type	Expiration	Vendor	Vendor Lot #	Vendor Cat Lot #	Volume Used	Volume Units
13C3 HFPO-DA_00004	13C3 HFPO-DA I.S. for HFPO	ASTD	08/28/18	Wellington Laboratories	M3HFPOADA0616M3HFPO-DA	1.00000		mL

Ok *Handwritten signature*  
 8/29/17


Dataset: Untitled

Last Altered: Tuesday, August 29, 2017 10:47:21 Mountain Daylight Time  
Printed: Tuesday, August 29, 2017 10:47:53 Mountain Daylight Time

Method: C:\MassLynx\8321.PRO\MethDB\hfpo.mdb 23 Aug 2017 10:19:52  
Library: C:\MassLynx\8321.PRO\CurveDB\hfpo17d24.cdb 24 Apr 2017 13:20:17

Sample Name: hfpo717H23083

HFPO IS 00004 MRM of 2 channels, ES-  
328.8 > 284.8  
2.895e+005

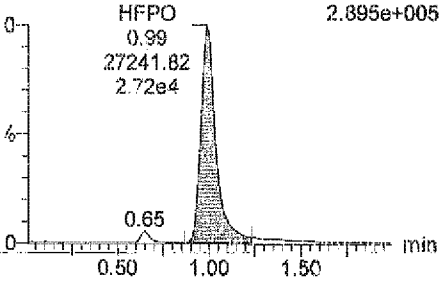


#	Name	Type	Std. Conc	RT	Area	IS Area	Response	Primar	ppb	%Dev
1	hfpo717H23083		10.000	0.99	27241.822	28795.438	0.946	bd	10.0	-0.4



Last Altered:   Tuesday, August 29, 2017 10:47:21 Mountain Daylight Time  
Created:        Tuesday, August 29, 2017 10:47:53 Mountain Daylight Time

Sample Name: hfpo717H23083

HFPO IS 00004     MRM of 2 channels, ES-  
                    331.8 > 286.8  
                    2.939e+005

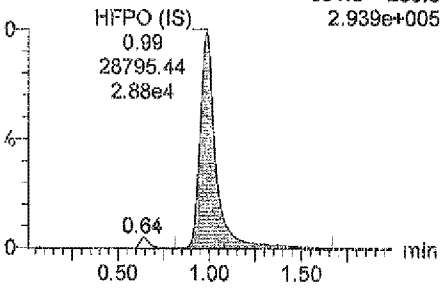

#	Name	Type	Std. Conc	RT	Area	IS Area	Response	Prima	ppb	%Dev
1	hfpo717H23083		1.000	0.99	28795.438		28795.438	bb	1.2	23.6

Reagent

---

**HFPO I.S.\_00007**

## Preliminary Report

TestAmerica Denver  
Internal Standard Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171212-65681.b\hfpo717L12074.d  
Lims ID: HFPO IS 00007  
Client ID:  
Sample Type: CCV  
Inject. Date: 12-Dec-2017 15:02:32 ALS Bottle#: 25 Worklist Smp#: 74  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: HFPO IS 00007  
Misc. Info.: HFPO17L12  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Sublist: chrom-HFPO\*sub1  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171212-65681.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 12-Dec-2017 15:48:38 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1: Det: F1:MRM  
Process Host: XAWRK024

## Averaged ICal Samples:

\\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10026.d  
\\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10027.d  
\\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10028.d  
\\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10029.d  
\\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10030.d  
\\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10031.d  
\\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10032.d  
\\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d

## Area Recoveries, Detector: F1:MRM

Compound	Average Standard	Lower Limit	Upper Limit	Sample	% Rec
* 2 13C3 HFPO-DA (IS)	731446	365723	1462892	740105	101.18

## RT Recoveries

Compound	Average Standard	Lower Limit	Upper Limit	Sample	DLT(min.)	% Diff
* 2 13C3 HFPO-DA (IS)	0.880	0.380	1.380	1.056	-0.176	19.997

AREA UPPER LIMIT = +100% of internal standard area.

AREA LOWER LIMIT = - 50% of internal standard area.

RT UPPER LIMIT = + 0.500 minutes of internal standard RT.

RT LOWER LIMIT = - 0.500 minutes of internal standard RT.

Reagent

---

**HFPO-DA\_00003**



# WELLINGTON LABORATORIES

## CERTIFICATE OF ANALYSIS DOCUMENTATION

**PRODUCT CODE:**

HFPO-DA

**LOT NUMBER:** HFPODA0213

**COMPOUND:**

2,5,3,3-Tetrafluoro-2-(1,1,2,2,3,3,3-heptafluoropropoxy)-propanoic acid

**STRUCTURE:**

**CAS #:**

13252-13-6



**MOLECULAR FORMULA:**

C<sub>8</sub>H<sub>5</sub>F<sub>11</sub>O<sub>5</sub>

**MOLECULAR WEIGHT:**

330.05

**CONCENTRATION:**

50 ± 2.5 µg/ml

**SOLVENT(S):**

Methanol

**CHEMICAL PURITY:**

>98%

**LAST TESTED:** (mm/dd/yyyy)

02/05/2014

**EXPIRY DATE:** (mm/dd/yyyy)

Stability studies ongoing

**RECOMMENDED STORAGE:**

Store ampoule in a cool, dark place

**DOCUMENTATION DATA ATTACHED:**

Figure 1: LC/MS Data (TIC and Mass Spectrum)

Figure 2: LC/MS/MS Data (Selected MRM Transitions)

**ADDITIONAL INFORMATION:**

See page 2 for further details.

**FOR LABORATORY USE ONLY: NOT FOR HUMAN OR DRUG USE**

**Certified By:**

B.G. Chittim

**Date:**

02/13/2014  
(mm/dd/yyyy)

Wellington Laboratories Inc., 345 Southgate Dr. Guelph ON N1G 3M5 CANADA  
519-822-2436 • Fax: 519-822-2849 • Info@well-labs.com

#### INTENDED USE:

The products prepared by Wellington Laboratories Inc. are for laboratory use only. They are designed to be used as reference standards for the identification and/or quantification of specific chemical compound(s).

#### HAZARDS:

This product should only be used by qualified personnel familiar with its potential hazards and trained in the handling of hazardous chemicals. Due care should be exercised to prevent unnecessary human contact or ingestion. All procedures should be carried out in a well-functioning fume hood and suitable gloves, eye protection and clothing should be worn at all times. Waste should be disposed of according to national and regional regulations. Material Safety Data Sheets (MSDSs) are available upon request.

#### SYNTHESIS / CHARACTERIZATION:

Where possible, all of our products are synthesized using single-product, unambiguous routes. They are then characterized, and their structures and purities confirmed, using a combination of the most relevant techniques, such as NMR, GC/MS, LC/MS/MS, x-ray crystallography and melting point. Isotopic purities of mass-labeled compounds are also confirmed using HRGC/HRMS and/or LC/MS/MS.

#### HOMOGENEITY:

Prior to solution preparation, crystalline material is tested for homogeneity using a variety of techniques (as stated above) and its solubility in a given solvent is taken into consideration. Duplicate solutions of a new product are prepared from the same crystalline lot and, after the addition of an appropriate internal standard, they are compared by GC/MS and/or LC/MS/MS. The relative response factors of the analyte of interest in each solution are required to be <5% RSD. New solution lots of existing products are compared to older lots in the same manner, which further confirms the homogeneity of the crystalline material as well as the stability and homogeneity of the solutions in the storage containers.

#### UNCERTAINTY:

The maximum combined relative standard uncertainty of our reference standard solutions is calculated using the following equation:

The combined relative standard uncertainty,  $u_c(y)$ , of a value  $y$  and the uncertainty of the independent parameters  $x_1, x_2, \dots, x_n$  on which it depends is:

$$u_c(y(x_1, x_2, \dots, x_n)) = \sqrt{\sum_{i=1}^n \left( \frac{\partial y}{\partial x_i} u(x_i) \right)^2}$$

where  $u$  is expressed as a relative standard uncertainty of the individual parameter.

The individual uncertainties taken into account include those associated with weights (calibration of the balance) and volumes (calibration of the volumetric glassware). An expanded maximum combined percent relative uncertainty of  $\pm 5\%$  (calculated with a coverage factor of 2 and a level of confidence of 95%) is stated on the Certificate of Analysis for all our products.

#### TRACEABILITY:

All reference standard solutions are traceable to specific crystalline lots. The microbalances used for solution preparation are regularly tested by an external, ISO/IEC 17025:2005 accredited calibration company. In addition, their calibration is verified prior to each weighing using NIST and/or NRC traceable external weights. All volumetric glassware used is of Class A tolerance and has been tested according to the appropriate ASTM procedures, which are ultimately traceable to NIST. For certain products, traceability to international interlaboratory studies has also been established.

#### EXPIRY DATE / PERIOD OF VALIDITY:


Ongoing stability studies of this product have demonstrated stability in its composition and concentration for the period of time specified by the expiry date in the unopened ampoule. Monitoring for any degradation or change in concentration of the listed analyte(s) is performed on a routine basis.

#### LIMITED WARRANTY:

At the time of shipment, all products are warranted to be free of defects in material and workmanship and to conform to the stated technical and purity specifications.



#### QUALITY MANAGEMENT:

This product was produced using a Quality Management System registered to ISO 9001:2008 by SAI Global, ISO/IEC 17025:2005 by the Canadian Association for Laboratory Accreditation Inc. (CALA, A 1226), and ISO GUIDE 34:2009 by ACLASS (certificate number AB-1522).

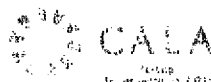

\*For additional information or assistance concerning this or any other products from Wellington Laboratories Inc., please visit our website at [www.well-labs.com](http://www.well-labs.com) or contact us directly at [info@well-labs.com](mailto:info@well-labs.com)\*\*

**Figure 1: HFPO-DA; LC/MS Data (TIC and Mass Spectrum)**

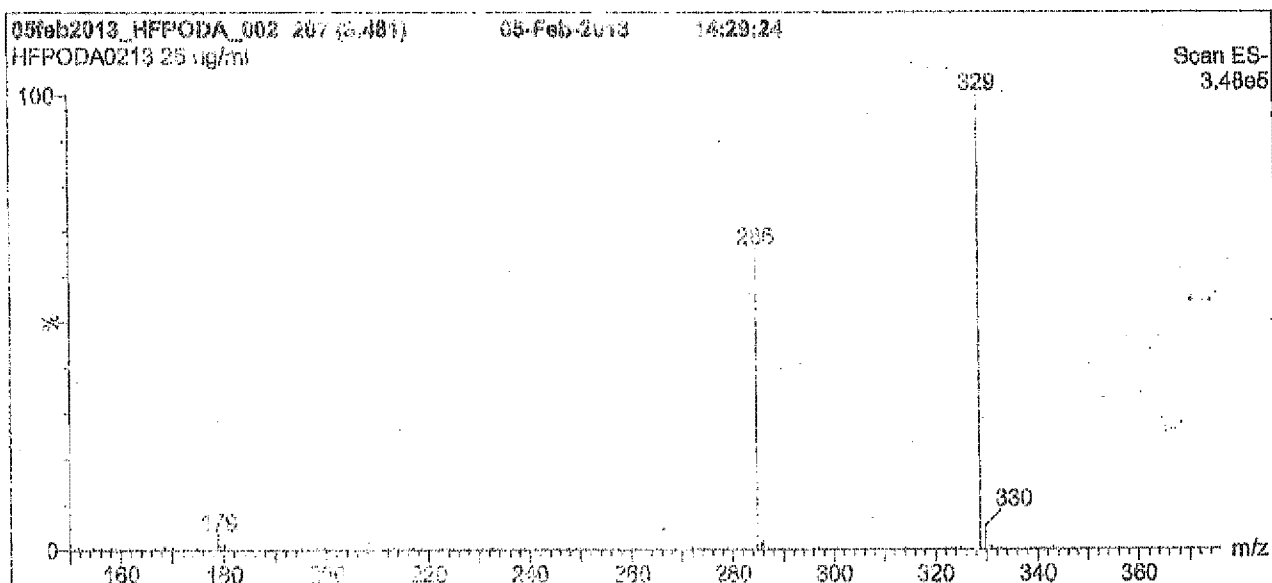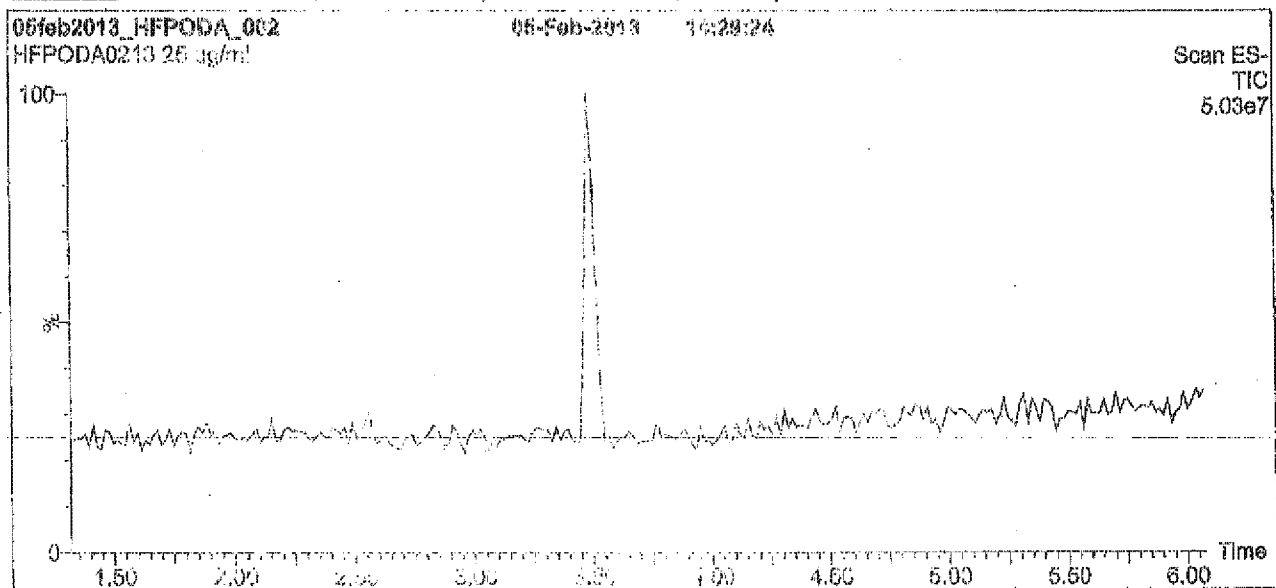

**Conditions for Figure 1:**

**LC:** Waters Acquity UPLC Performance LC  
**MS:** Micromass Quattro micro API MS

**Chromatographic Conditions:**

**Column:** Kinetex PFP  
2.6 µm, 4.6 x 100 mm

**Mobile phase:** Gradient  
Start: 40% (00:20) 0.1% H<sub>3</sub>PO<sub>4</sub> / 60% H<sub>2</sub>O  
(both in 10 mM NH<sub>4</sub> OAc buffer)  
Ramp to 80% organic over 9 min and hold for 1 min  
before returning to initial conditions in 0.5 min.  
Time: 11 min


**Flow:** 800 µl/min

**MS Parameters:**

**Experiment:** Full Scan (150 - 850 amu)

**Source:** Electrospray (negative)  
**Capillary Voltage (kV):** 3.00  
**Cone Voltage (V):** 9.00  
**Cone Gas Flow (l/hr):** 50  
**Desolvation Gas Flow (l/hr):** 750

**Figure 2: HFPO-DA: LC/MS/MS Data (Selected MRM Transitions)**

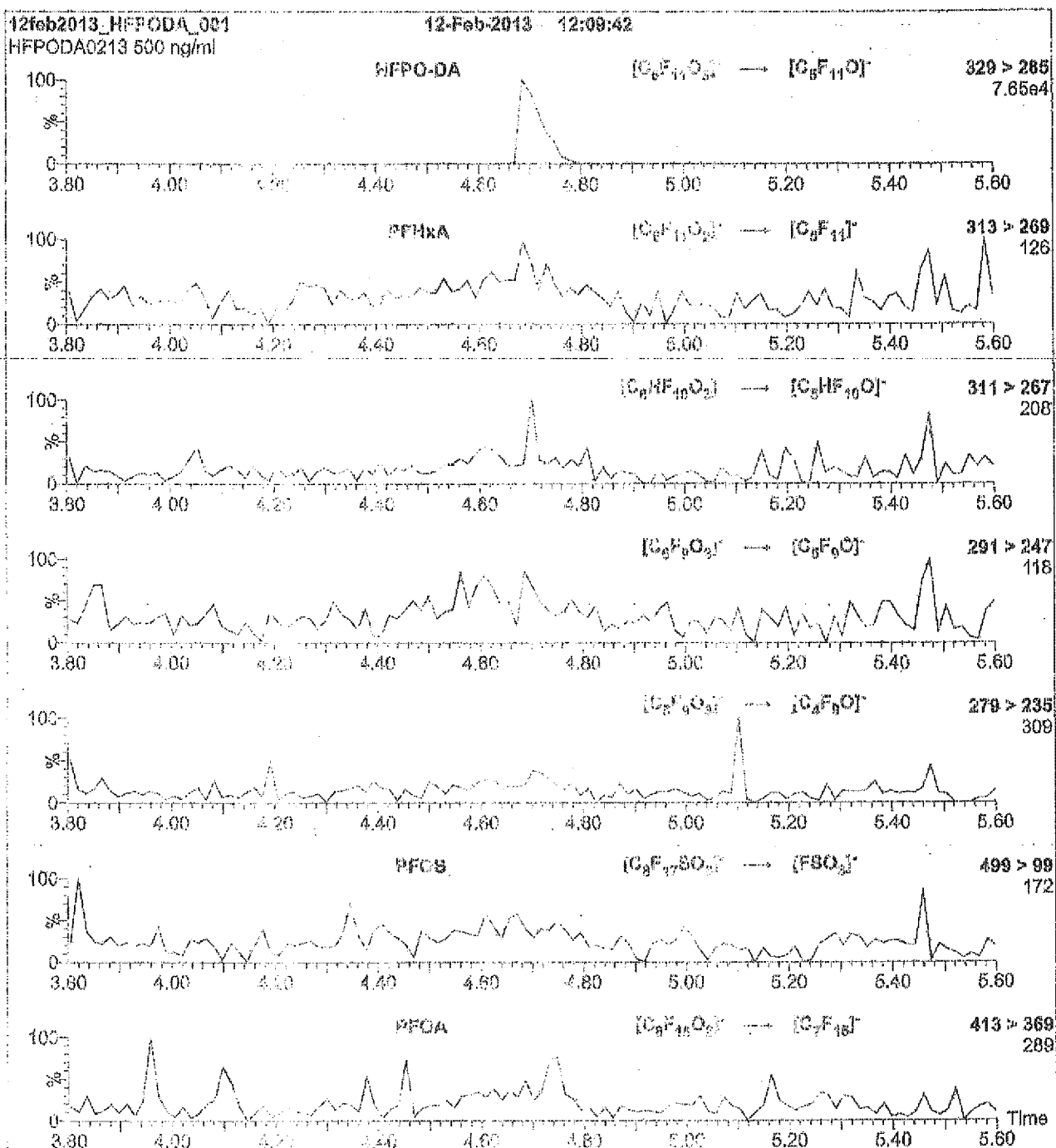

**Conditions for Figure 2:**

**Injection:** Direct loop injection  
10  $\mu$ l (500 ng/ml HFPO-DA)

**Mobile phase:** Isocratic 80% (80:20 MeOH:ACN) / 20%  $H_2O$   
(both with 10 mM NH<sub>4</sub>OAc buffer)

**Flow:** 300  $\mu$ l/min

**MS Parameters:**

Collision Gas (mbar) = 3.87e-3  
Collision Energy (eV) = 5



Reagent

---

**HFPO-DA\_00004**



# WELLINGTON LABORATORIES

## CERTIFICATE OF ANALYSIS DOCUMENTATION

**PRODUCT CODE:**

HFPO-DA

**LOT NUMBER:**


HFPODA0717

**COMPOUND:**

2,3,3,3-Tetrafluoro-2-(1,1,2,2,3,3,3-heptafluoropropoxy)-propanoic acid

**STRUCTURE:****CAS #:**

13252-13-6

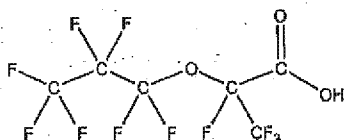**MOLECULAR FORMULA:** $C_6H_5F_{11}O_3$ **MOLECULAR WEIGHT:**

330.05

**CONCENTRATION:** $50 \pm 2.5 \mu\text{g/ml}$ **SOLVENT(S):**

Methanol

**CHEMICAL PURITY:**

&gt;98%

**LAST TESTED:** (mm/dd/yyyy)

07/13/2017

**EXPIRY DATE:** (mm/dd/yyyy)

07/13/2020

**RECOMMENDED STORAGE:**

Store ampoule in a cool, dark place

**DOCUMENTATION/ DATA ATTACHED:**

Figure 1: LC/MS Data (TIC and Mass Spectrum)

Figure 2: LC/MS/MS Data (Selected MRM Transitions)

**ADDITIONAL INFORMATION:**

- See page 2 for further details.
- Product is commercially known as GenX.

**FOR LABORATORY USE ONLY: NOT FOR HUMAN OR DRUG USE**

Certified By:

  
B.G. Chittim, General Manager

Date:

07/14/2017  
(mm/dd/yyyy)

Wellington Laboratories Inc., 345 Southgate Dr. Guelph ON N1G 3M5 CANADA  
519-822-2436 • Fax: 519-822-2849 • [info@well-labs.com](mailto:info@well-labs.com)

**INTENDED USE:**

The products prepared by Wellington Laboratories Inc. are for laboratory use only. This certified reference material (CRM) was designed to be used as a standard for the identification and/or quantification of the specific chemical compound it contains.

**HAZARDS:**

This product should only be used by qualified personnel familiar with its potential hazards and trained in the handling of hazardous chemicals. Due care should be exercised to prevent unnecessary human contact or ingestion. All procedures should be carried out in a well-functioning fume hood and suitable gloves, eye protection, and clothing should be worn at all times. Waste should be disposed of according to national and regional regulations. Safety Data Sheets (SDSs) are available upon request.

**SYNTHESIS / CHARACTERIZATION:**

Where possible, all of our products are synthesized using single-product unambiguous routes. They are then characterized, and their structures and purities confirmed, using a combination of the most relevant techniques, such as NMR, GC/MS, LC/MS/MS, SFC/UV/MS/MS, x-ray crystallography, and melting point. Isotopic purities of mass-labelled compounds are also confirmed using HRGC/HRMS and/or LC/MS/MS.

**HOMOGENEITY:**

Prior to solution preparation, crystalline material is tested for homogeneity using a variety of techniques (as stated above) and its solubility in a given diluent is taken into consideration. Duplicate solutions of a new product are prepared from the same crystalline lot and, after the addition of an appropriate internal standard, they are compared by GC/MS, LC/MS/MS and/or SFC/UV/MS/MS. The relative response factors of the analyte of interest in each solution are required to be <5% RSD. New solution lots of existing products are compared to older lots in the same manner, which further confirms the homogeneity of the crystalline material as well as the stability and homogeneity of the solutions in the storage containers. In order to maintain the integrity of the assigned value(s), and associated uncertainty, the dilution or injection of a subsample of this product should be performed using calibrated measuring equipment.

**UNCERTAINTY:**

The maximum combined relative standard uncertainty of our reference standard solutions is calculated using the following equation:

The combined relative standard uncertainty,  $u_c(y)$ , of a value  $y$  and the uncertainty of the independent parameters

$x_1, x_2, \dots, x_n$  on which it depends is:

$$u_c(y(x_1, x_2, \dots, x_n)) = \sqrt{\sum_{i=1}^n u(y, x_i)^2}$$

where  $x$  is expressed as a relative standard uncertainty of the individual parameter.

The individual uncertainties taken into account include those associated with weights (calibration of the balance) and volumes (calibration of the volumetric glassware). An expanded maximum combined percent relative uncertainty of  $\pm 5\%$  (calculated with a coverage factor of 2 and a level of confidence of 95%) is stated on the Certificate of Analysis for all of our products.

**TRACEABILITY:**

All reference standard solutions are traceable to specific crystalline lots. The microbalances used for solution preparation are regularly tested by an external ISO/IEC 17025 accredited calibration company. In addition, their calibration is verified prior to each weighing using calibrated NIST and/or NRC traceable external weights. All volumetric glassware used is calibrated, of Class A tolerance, and has been tested according to the appropriate ASTM procedures, which are ultimately traceable to NIST. For certain products, traceability to International Interlaboratory studies has also been established.

**EXPIRY DATE / PERIOD OF VALIDITY:**

Ongoing stability studies of this product have demonstrated stability in its composition and concentration, until the specified expiry date, in the unopened ampoule. Monitoring for any degradation or change in concentration of the listed analyte(s) is performed on a routine basis.

**LIMITED WARRANTY:**

At the time of shipment, all products are warranted to be free of defects in material and workmanship and to conform to the stated technical and purity specifications.



**QUALITY MANAGEMENT:**

This product was produced using a Quality Management System registered to the latest versions of ISO 9001 by SAI Global, ISO/IEC 17025 by the Canadian Association for Laboratory Accreditation Inc. (CALA; A 1226), and ISO GUIDE 34 by ANSI-ASQ National Accreditation Board (ANAB; AR-1523).



\*\*For additional information or assistance concerning this or any other products from Wellington Laboratories Inc., please visit our website at [www.well-labs.com](http://www.well-labs.com) or contact us directly at [info@well-labs.com](mailto:info@well-labs.com)\*\*

**Figure 1: HFPO-DA; LC/MS Data (TIC and Mass Spectrum)**

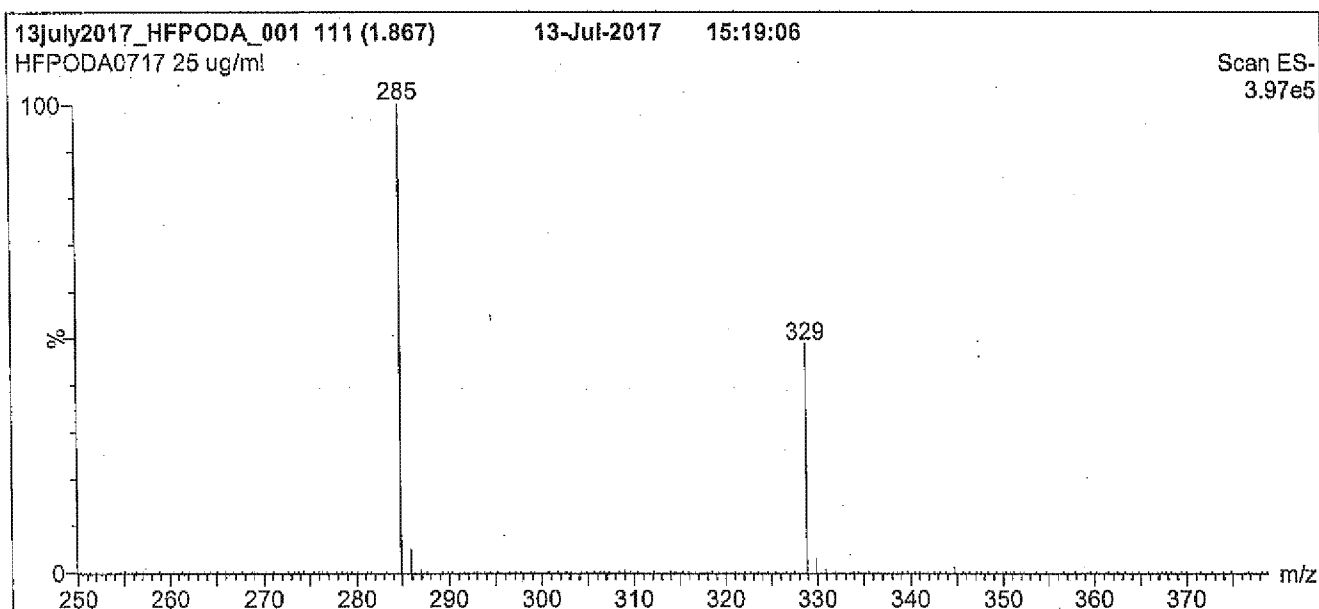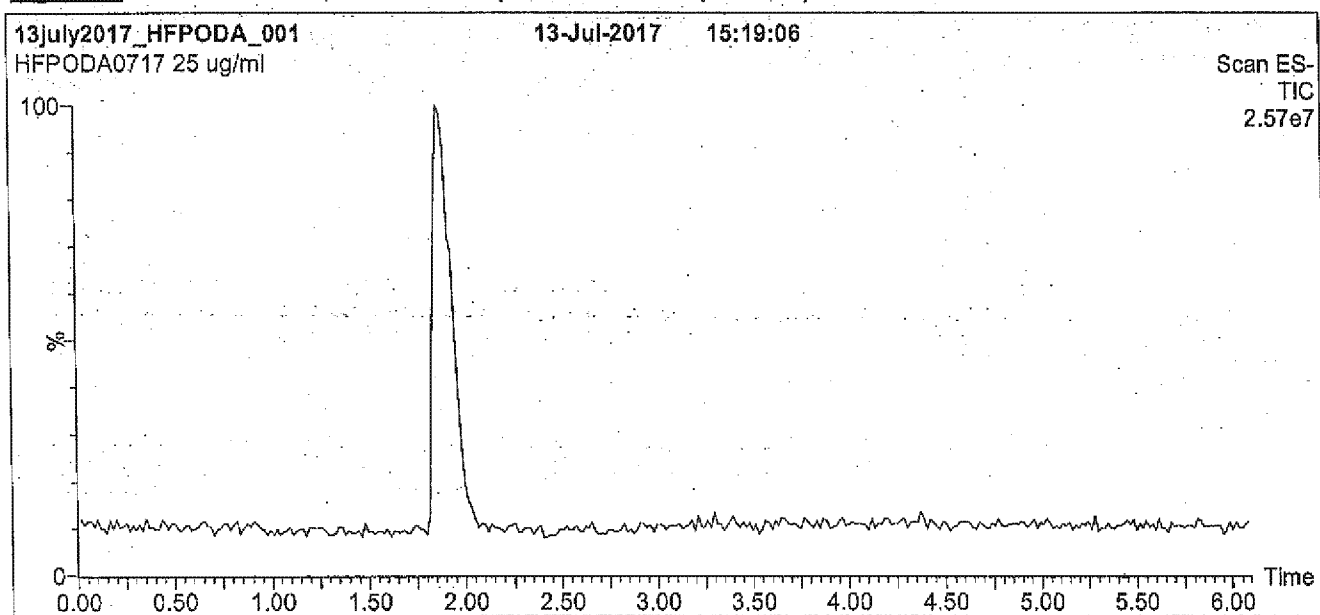

**Conditions for Figure 1:**

**LC:** Waters Acquity Ultra Performance LC  
**MS:** Micromass Quattro *micro* API MS

**Chromatographic Conditions**

Column: Acquity UPLC BEH Shield RP<sub>18</sub>  
1.7  $\mu$ m, 2.1 x 100 mm

Mobile phase: Gradient

Start: 55% MeOH / 45% H<sub>2</sub>O with 10 mM NH<sub>4</sub>OAc buffer  
Ramp to 90% organic over 7.5 min and hold for 1.5 min  
before returning to initial conditions in 0.5 min.

Time: 10 min

Flow: 300  $\mu$ l/min

**MS Parameters**

Experiment: Full Scan (250 - 850 amu)

Source: Electrospray (negative)


Capillary Voltage (kV) = 3.00

Cone Voltage (V) = 10.00

Cone Gas Flow (l/hr) = 100

Desolvation Gas Flow (l/hr) = 700

**Figure 2: HFPO-DA; LC/MS/MS Data (Selected MRM Transitions)**

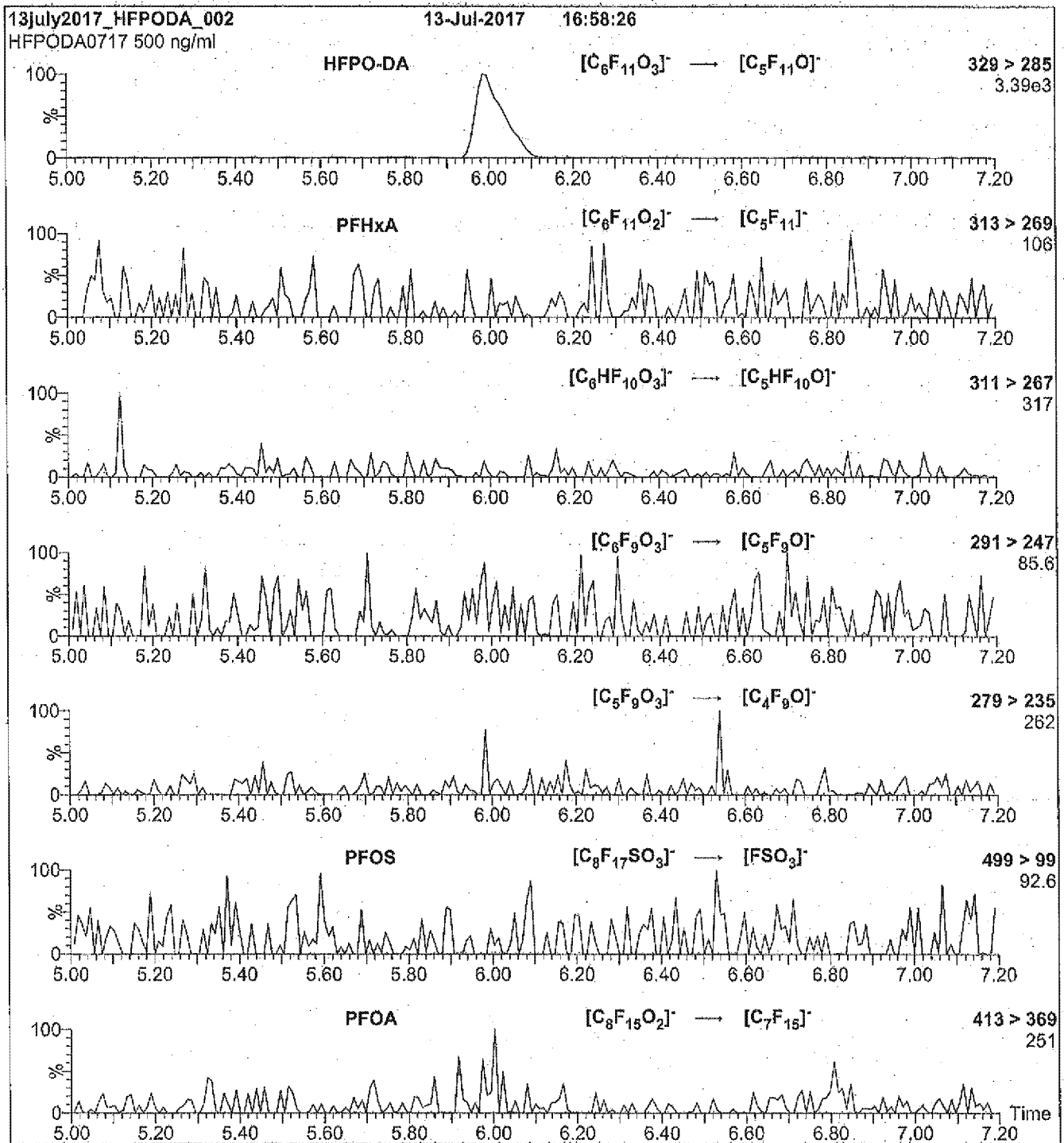

**Conditions for Figure 2:**

Injection: Direct loop Injection  
10  $\mu$ l (500 ng/ml HFPO-DA)

Mobile phase: Isocratic 80% MeOH / 20% H<sub>2</sub>O with 10 mM NH<sub>4</sub>OAc buffer

Flow: 300  $\mu$ l/min

**MS Parameters**

Collision Gas (mbar) = 3.20e-3  
Collision Energy (eV) = 5

# 8321A\_HFPO\_Du

---

HFPO-DA

FORM II  
LCMS SURROGATE RECOVERY

Lab Name: TestAmerica Denver

Job No.: 280-105698-1

SDG No.: \_\_\_\_\_

Matrix: Water

Level: Low

GC Column (1): Synergi Hyd ID: \_\_\_\_\_

Client Sample ID	Lab Sample ID	HFPODA #
FAY-D-6394CHKFT-W1-1-012218	280-105698-1	89
FAY-D-6394CHKFT-W1-1-012218-D	280-105698-2	88
FAY-D-6246CHKFT-W1-1-012218	280-105698-3	89
FAY-D-318BOONE-W1-1-012218	280-105698-4	92
FAY-D-41BOONE-W1-1-012218	280-105698-5	88
FAY-D-FB-012218	280-105698-6	92
FAY-D-7145BUTLE-W1-1-012218	280-105698-7	84
FAY-D-1515SCLLY-W1-1-012218	280-105698-8	88
FAY-D-7396SALIE-W1-1-012218	280-105698-9	86
FAY-D-7012NC87H-W1-1-012218	280-105698-10	86
	MB 280-402648/1-A	93
	LCS 280-402648/2-A	91
	LCSD 280-402648/3-A	90
	LLCS 280-402648/4-A	94
FAY-D-6394CHKFT-W1-1-012218 MS	280-105698-1 MS	86
FAY-D-6394CHKFT-W1-1-012218 DU	280-105698-1 DU	85
	DLCK 280-390728/12	102

HFPODA = 13C3 HFPO-DA

QC LIMITS  
50-200

# Column to be used to flag recovery values

FORM II 8321A

FORM III  
LCMS LAB CONTROL SAMPLE RECOVERY

Lab Name: TestAmerica Denver Job No.: 280-105698-1  
SDG No.: \_\_\_\_\_  
Matrix: Water Level: Low Lab File ID: hfpo718A25008.d  
Lab ID: LCS 280-402648/2-A Client ID: \_\_\_\_\_

COMPOUND	SPIKE ADDED (ug/L)	LCS CONCENTRATION (ug/L)	LCS % REC	QC LIMITS REC	#
HFPO-DA	0.200	0.193	96	70-130	

# Column to be used to flag recovery and RPD values



FORM III  
LCMS LAB CONTROL SAMPLE DUPLICATE RECOVERY

Lab Name: TestAmerica Denver Job No.: 280-105698-1  
SDG No.: \_\_\_\_\_  
Matrix: Water Level: Low Lab File ID: hfpo718A25009.d  
Lab ID: LCSD 280-402648/3-A Client ID: \_\_\_\_\_

COMPOUND	SPIKE ADDED (ug/L)	LCSD CONCENTRATION (ug/L)	LCSD % REC	% RPD	QC LIMITS		#
					RPD	REC	
HFPO-DA	0.200	0.190	95	1	20	70-130	

# Column to be used to flag recovery and RPD values

FORM III  
LCMS LOW LEVEL CONTROL SAMPLE RECOVERY

Lab Name: TestAmerica Denver Job No.: 280-105698-1  
SDG No.: \_\_\_\_\_  
Matrix: Water Level: Low Lab File ID: hfpo718A25010.d  
Lab ID: LLCS 280-402648/4-A Client ID: \_\_\_\_\_

COMPOUND	SPIKE ADDED (ug/L)	LLCS CONCENTRATION (ug/L)	LLCS % REC	QC LIMITS REC	#
HFPO-DA	0.0200	0.0173	87	70-130	

# Column to be used to flag recovery and RPD values

FORM III  
LCMS MATRIX SPIKE RECOVERY

Lab Name: TestAmerica Denver Job No.: 280-105698-1  
SDG No.: \_\_\_\_\_  
Matrix: Water Level: Low Lab File ID: hfpo718A25013.d  
Lab ID: 280-105698-1 MS Client ID: FAY-D-6394CHKFT-W1-1-012218 MS

COMPOUND	SPIKE ADDED (ug/L)	SAMPLE CONCENTRATION (ug/L)	MS CONCENTRATION (ug/L)	MS % REC	QC LIMITS REC	#
HFPO-DA	0.198	0.033	0.220	94	70-130	

# Column to be used to flag recovery and RPD values

FORM III  
LCMS DETECTION LIMIT CHECK STANDARD RECOVERY

Lab Name: TestAmerica Denver Job No.: 280-105698-1  
SDG No.: \_\_\_\_\_  
Matrix: Water Level: Low Lab File ID: hfpo717J10035.d  
Lab ID: DLCK 280-390728/12 Client ID: \_\_\_\_\_

COMPOUND	SPIKE ADDED (ug/L)	DLCK CONCENTRATION (ug/L)	DLCK % REC	QC LIMITS REC	#
HFPO-DA	0.250	<0.50	78	70-130	

# Column to be used to flag recovery and RPD values

FORM IV  
LCMS METHOD BLANK SUMMARY

Lab Name: <u>TestAmerica Denver</u>	Job No.: <u>280-105698-1</u>
SDG No.: _____	
Lab File ID: <u>hfpo718A25007.d</u>	Lab Sample ID: <u>MB 280-402648/1-A</u>
Matrix: <u>Water</u>	Date Extracted: <u>01/24/2018 15:20</u>
Instrument ID: <u>LC_LCMS7</u>	Date Analyzed: <u>01/25/2018 10:41</u>
Level: (Low/Med) <u>Low</u>	

THIS METHOD BLANK APPLIES TO THE FOLLOWING SAMPLES:

CLIENT SAMPLE ID	LAB SAMPLE ID	LAB FILE ID	DATE ANALYZED
	LCS 280-402648/2-A	hfpo718A25008.d	01/25/2018 10:44
	LCSD 280-402648/3-A	hfpo718A25009.d	01/25/2018 10:48
	LLCS 280-402648/4-A	hfpo718A25010.d	01/25/2018 10:51
FAY-D-6394CHKFT-W1-1-012218	280-105698-1	hfpo718A25011.d	01/25/2018 10:54
FAY-D-6394CHKFT-W1-1-012218 DU	280-105698-1 DU	hfpo718A25012.d	01/25/2018 10:57
FAY-D-6394CHKFT-W1-1-012218 MS	280-105698-1 MS	hfpo718A25013.d	01/25/2018 11:01
FAY-D-6394CHKFT-W1-1-012218-D	280-105698-2	hfpo718A25015.d	01/25/2018 11:07
FAY-D-6246CHKFT-W1-1-012218	280-105698-3	hfpo718A25016.d	01/25/2018 11:10
FAY-D-318BOONE-W1-1-012218	280-105698-4	hfpo718A25017.d	01/25/2018 11:14
FAY-D-41BOONE-W1-1-012218	280-105698-5	hfpo718A25018.d	01/25/2018 11:17
FAY-D-FB-012218	280-105698-6	hfpo718A25019.d	01/25/2018 11:20
FAY-D-7145BUTLE-W1-1-012218	280-105698-7	hfpo718A25020.d	01/25/2018 11:23
FAY-D-1515SCLLY-W1-1-012218	280-105698-8	hfpo718A25021.d	01/25/2018 11:27
FAY-D-7396SALIE-W1-1-012218	280-105698-9	hfpo718A25023.d	01/25/2018 11:33
FAY-D-7012NC87H-W1-1-012218	280-105698-10	hfpo718A25024.d	01/25/2018 11:36



TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25011.d  
Lims ID: 280-105698-D-1-A  
Client ID: FAY-D-6394CHKFT-W1-1-012218  
Sample Type: Client  
Inject. Date: 25-Jan-2018 10:54:35 ALS Bottle#: 15 Worklist Smp#: 8  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: 280-105698-D-1-A  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:16 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyera

Date: 25-Jan-2018 15:25:52

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 &gt; 286.8 0.961 0.961 0.0 1.000 648766 8.87 1292

\* 2 13C3 HFPO-DA (IS)

331.8 &gt; 286.8 0.961 0.961 0.0 648766 10.0 1292

1 Perfluoro(2-propoxypropanoic) acid

328.8 &gt; 284.8 0.961 0.988 -0.027 1.000 130017 1.76 44.7 M

## QC Flag Legend

Review Flags

M - Manually Integrated

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25011.d

Injection Date: 25-Jan-2018 10:54:35

Instrument ID: LC\_LCMS7

Lims ID: 280-105698-D-1-A

Lab Sample ID: 280-105698-1

Client ID: FAY-D-6394CHKFT-W1-1-012218

Operator ID: JBH

ALS Bottle#: 15

Worklist Smp#: 8

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid (M)

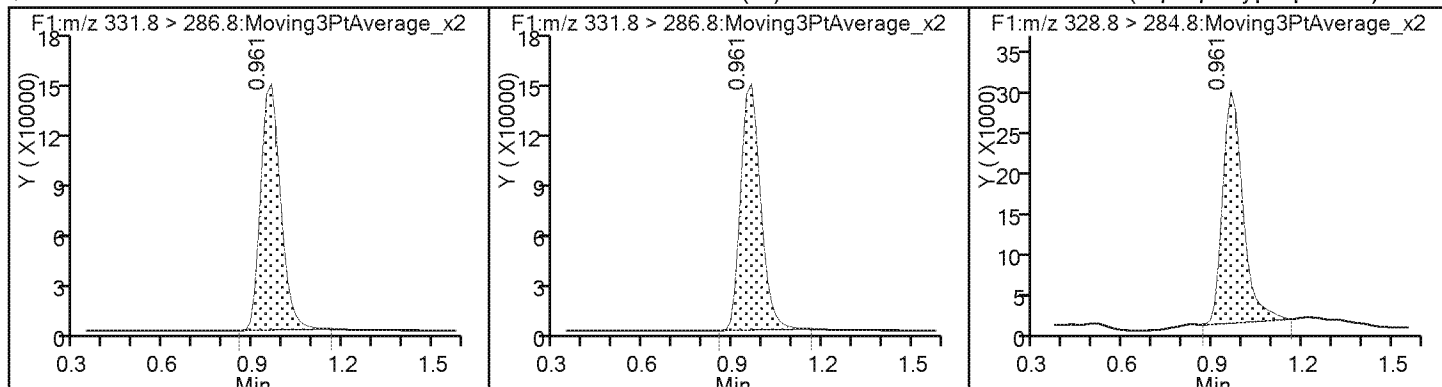



TestAmerica Denver  
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25011.d  
Lims ID: 280-105698-D-1-A  
Client ID: FAY-D-6394CHKFT-W1-1-012218  
Sample Type: Client  
Inject. Date: 25-Jan-2018 10:54:35 ALS Bottle#: 15 Worklist Smp#: 8  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: 280-105698-D-1-A  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:16 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyera

Date: 25-Jan-2018 15:25:52

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	8.87	88.70

## TestAmerica Denver


Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25011.d  
Injection Date: 25-Jan-2018 10:54:35 Instrument ID: LC\_LCMS7  
Lims ID: 280-105698-D-1-A Lab Sample ID: 280-105698-1  
Client ID: FAY-D-6394CHKFT-W1-1-012218  
Operator ID: JBH ALS Bottle#: 15 Worklist Smp#: 8  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Method: HFPO Limit Group: LC - 8321A\_HFPO\_Du  
Column: Detector F1:MRM

## 1 Perfluoro(2-propoxypropanoic) acid, CAS: 13252-13-6

Signal: 1


RT: 0.96  
Area: 123636  
Amount: 1.667065  
Amount Units: ug/l

## Processing Integration Results

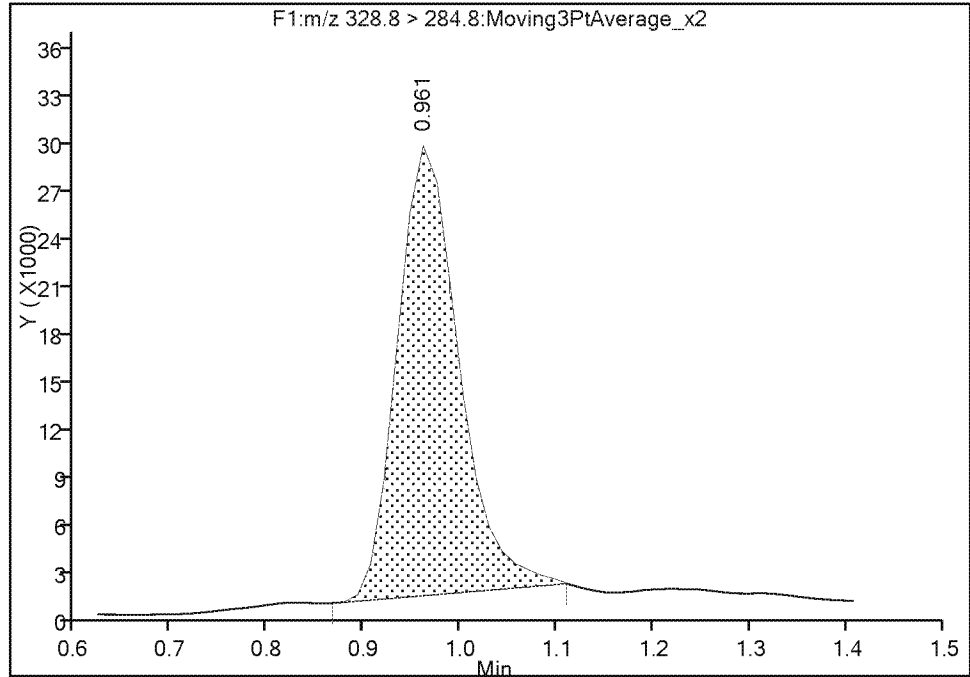

RT: 0.96  
Area: 130017  
Amount: 1.764246  
Amount Units: ug/l

## Manual Integration Results

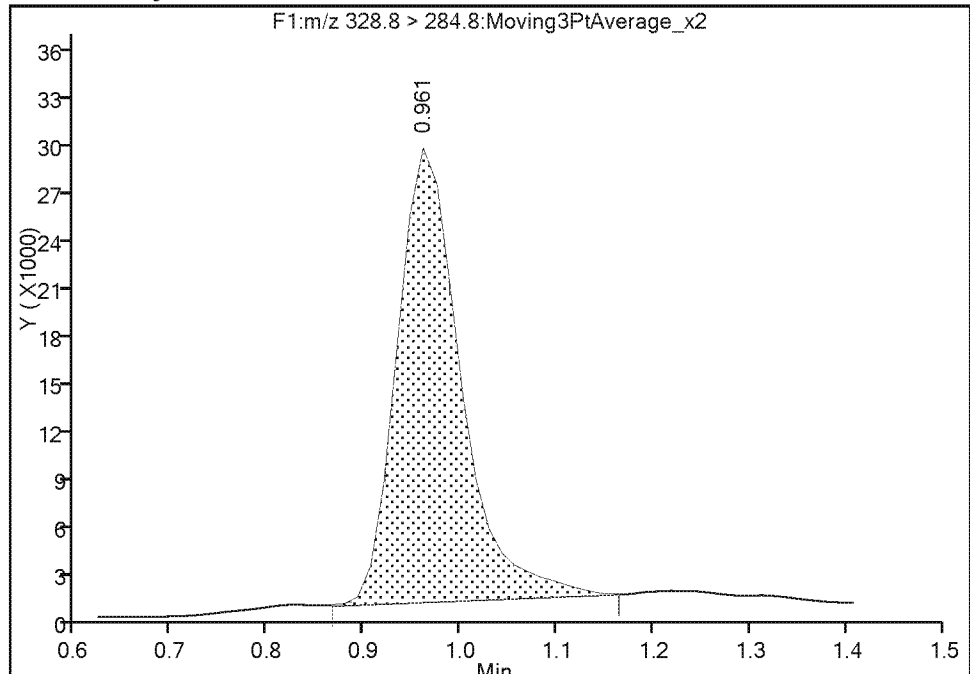

Reviewer: meyer, 25-Jan-2018 15:27:17

Audit Action: Manually Integrated

Audit Reason: Baseline

## FORM I

Lab Name: <u>TestAmerica Denver</u>	Job No.: <u>280-105698-1</u>
SDG No.: _____	
Client Sample ID: <u>FAY-D-6394CHKFT-W1-1-0122</u> <u>18-D</u>	Lab Sample ID: <u>280-105698-2</u>
Matrix: <u>Water</u>	Lab File ID: <u>hfpo718A25015.d</u>
Analysis Method: <u>8321A</u>	Date Collected: <u>01/22/2018</u> <u>16:19</u>
Extraction Method: <u>3535</u>	Date Extracted: <u>01/24/2018</u> <u>15:20</u>
Sample wt/vol: <u>250.7(mL)</u>	Date Analyzed: <u>01/25/2018</u> <u>11:07</u>
Con. Extract Vol.: <u>5(mL)</u>	Dilution Factor: <u>1</u>
Injection Volume: <u>20(uL)</u>	GC Column: <u>Synergi Hydro</u> ID: _____
% Moisture: _____	GPC Cleanup: (Y/N) <u>N</u>
Analysis Batch No.: <u>402806</u>	Units: <u>ug/L</u>

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	0.032		0.010	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	88		50-200

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25015.d  
Lims ID: 280-105698-B-2-A  
Client ID: FAY-D-6394CHKFT-W1-1-012218-D  
Sample Type: Client  
Inject. Date: 25-Jan-2018 11:07:33 ALS Bottle#: 18 Worklist Smp#: 12  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: 280-105698-B-2-A  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:23 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyera

Date: 25-Jan-2018 15:28:27

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 &gt; 286.8 1.002 0.961 0.041 1.000 640461 8.76 1719

\* 2 13C3 HFPO-DA (IS)

331.8 &gt; 286.8 1.002 0.961 0.041 640461 10.0 1719

1 Perfluoro(2-propoxypropanoic) acid

328.8 &gt; 284.8 1.015 0.988 0.027 1.000 117567 1.60 32.0 M

## QC Flag Legend

Review Flags

M - Manually Integrated

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25015.d

Injection Date: 25-Jan-2018 11:07:33

Instrument ID: LC\_LCMS7

Lims ID: 280-105698-B-2-A

Lab Sample ID: 280-105698-2

Client ID: FAY-D-6394CHKFT-W1-1-012218-D

Operator ID: JBH

ALS Bottle#: 18

Worklist Smp#: 12

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid (M)

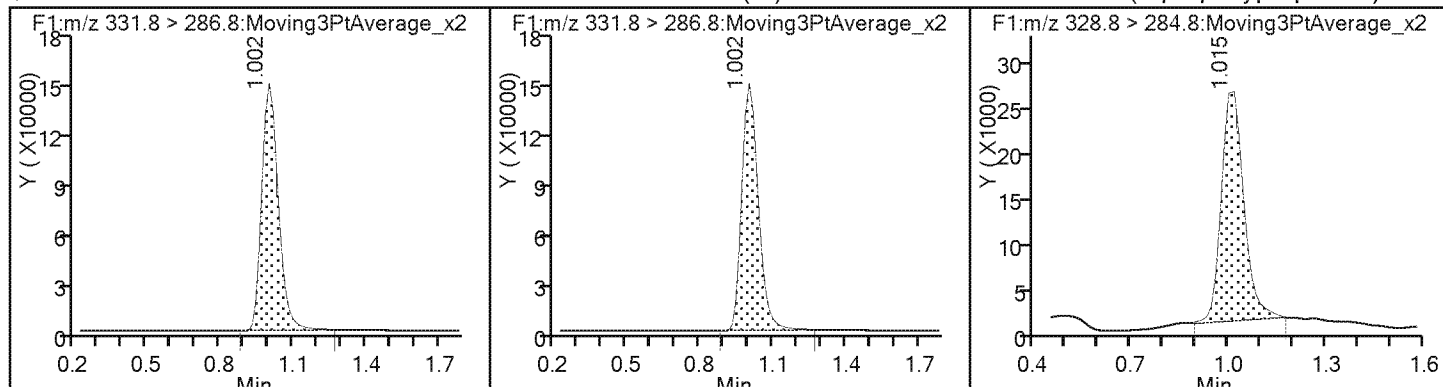

TestAmerica Denver  
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25015.d  
Lims ID: 280-105698-B-2-A  
Client ID: FAY-D-6394CHKFT-W1-1-012218-D  
Sample Type: Client  
Inject. Date: 25-Jan-2018 11:07:33 ALS Bottle#: 18 Worklist Smp#: 12  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: 280-105698-B-2-A  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:23 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyera

Date: 25-Jan-2018 15:28:27

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	8.76	87.56

## TestAmerica Denver


Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25015.d  
Injection Date: 25-Jan-2018 11:07:33 Instrument ID: LC\_LCMS7  
Lims ID: 280-105698-B-2-A Lab Sample ID: 280-105698-2  
Client ID: FAY-D-6394CHKFT-W1-1-012218-D  
Operator ID: JBH ALS Bottle#: 18 Worklist Smp#: 12  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Method: HFPO Limit Group: LC - 8321A\_HFPO\_Du  
Column: Detector F1:MRM

## 1 Perfluoro(2-propoxypropanoic) acid, CAS: 13252-13-6

Signal: 1


RT: 1.02  
Area: 112707  
Amount: 1.522878  
Amount Units: ug/l

## Processing Integration Results

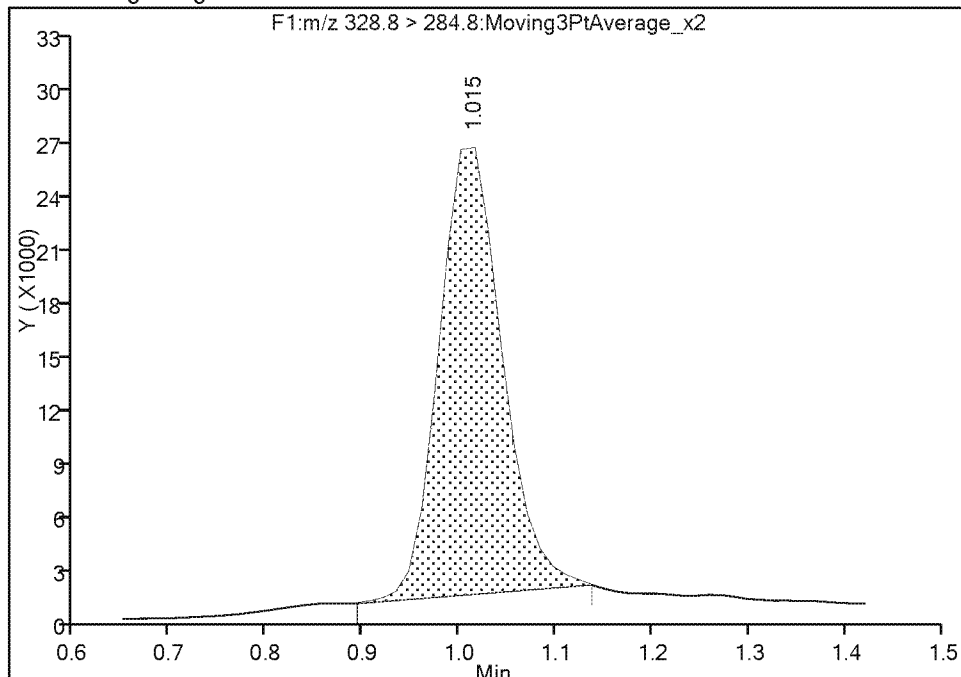

RT: 1.02  
Area: 117567  
Amount: 1.597854  
Amount Units: ug/l

## Manual Integration Results

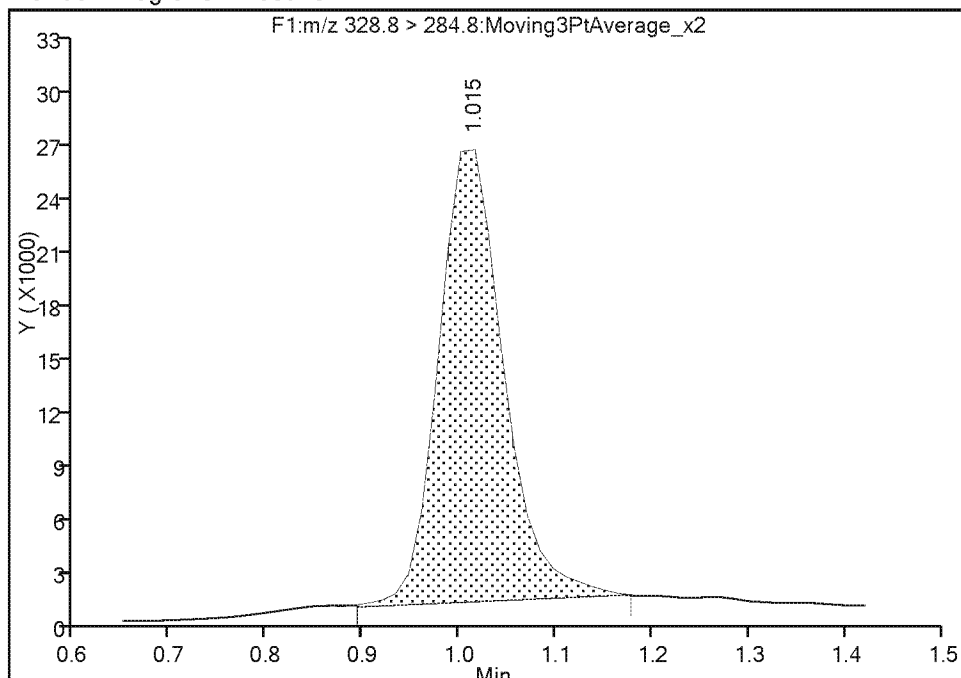

Reviewer: meyera, 25-Jan-2018 15:27:33

Audit Action: Manually Integrated

Audit Reason: Baseline

## FORM I

Lab Name: TestAmerica Denver Job No.: 280-105698-1

SDG No.:

Client Sample ID: FAY-D-6246CHKFT-W1-1-0122    Lab Sample ID: 280-105698-3

Matrix: Water Lab File ID: hfpo718A25016.d

Analysis Method: 8321A Date Collected: 01/22/2018 14:56

Extraction Method: 3535 Date Extracted: 01/24/2018 15:20

Sample wt/vol: 254.1(mL) Date Analyzed: 01/25/2018 11:10

Con. Extract Vol.: 5(mL) Dilution Factor: 1

Injection Volume: 20(uL) GC Column: Synergi Hydro ID:

% Moisture: GPC Cleanup: (Y/N) N

Analysis Batch No.: 402806 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	0.052		0.010	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	89		50-200



TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25016.d  
Lims ID: 280-105698-D-3-A  
Client ID: FAY-D-6246CHKFT-W1-1-012218  
Sample Type: Client  
Inject. Date: 25-Jan-2018 11:10:47 ALS Bottle#: 19 Worklist Smp#: 13  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: 280-105698-D-3-A  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:23 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyera

Date: 25-Jan-2018 15:28:32

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 &gt; 286.8 0.961 0.961 0.0 1.000 648364 8.86 1464

\* 2 13C3 HFPO-DA (IS)

331.8 &gt; 286.8 0.961 0.961 0.0 648364 10.0 1464

1 Perfluoro(2-propoxypropanoic) acid

328.8 &gt; 284.8 0.961 0.988 -0.027 1.000 187553 2.64 86.0

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25016.d

Injection Date: 25-Jan-2018 11:10:47

Instrument ID: LC\_LCMS7

Lims ID: 280-105698-D-3-A

Lab Sample ID: 280-105698-3

Client ID: FAY-D-6246CHKFT-W1-1-012218

Operator ID: JBH

ALS Bottle#: 19

Worklist Smp#: 13

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

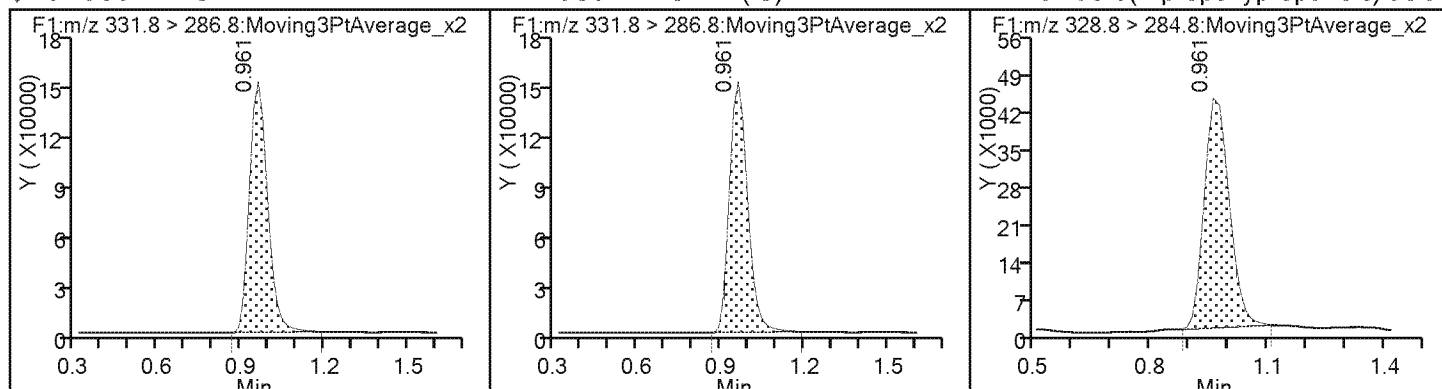

TestAmerica Denver  
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25016.d  
Lims ID: 280-105698-D-3-A  
Client ID: FAY-D-6246CHKFT-W1-1-012218  
Sample Type: Client  
Inject. Date: 25-Jan-2018 11:10:47 ALS Bottle#: 19 Worklist Smp#: 13  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: 280-105698-D-3-A  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:23 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyera

Date: 25-Jan-2018 15:28:32

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	8.86	88.64

Lab Name: <u>TestAmerica Denver</u>	Job No.: <u>280-105698-1</u>
SDG No.: _____	
Client Sample ID: <u>FAY-D-318BOONE-W1-1-01221</u> <u>8</u>	Lab Sample ID: <u>280-105698-4</u>
Matrix: <u>Water</u>	Lab File ID: <u>hfpo718A25017.d</u>
Analysis Method: <u>8321A</u>	Date Collected: <u>01/22/2018</u> <u>14:15</u>
Extraction Method: <u>3535</u>	Date Extracted: <u>01/24/2018</u> <u>15:20</u>
Sample wt/vol: <u>245.1(mL)</u>	Date Analyzed: <u>01/25/2018</u> <u>11:14</u>
Con. Extract Vol.: <u>5(mL)</u>	Dilution Factor: <u>1</u>
Injection Volume: <u>20(uL)</u>	GC Column: <u>Synergi Hydro</u> ID: _____
% Moisture: _____	GPC Cleanup: (Y/N) <u>N</u>
Analysis Batch No.: <u>402806</u>	Units: <u>ug/L</u>

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	0.044		0.010	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	92		50-200

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25017.d  
Lims ID: 280-105698-A-4-A  
Client ID: FAY-D-318BOONE-W1-1-012218  
Sample Type: Client  
Inject. Date: 25-Jan-2018 11:14:03 ALS Bottle#: 20 Worklist Smp#: 14  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: 280-105698-A-4-A  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:23 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyer

Date: 25-Jan-2018 15:28:34

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 &gt; 286.8 0.947 0.961 -0.014 1.000 671887 9.19 1479

\* 2 13C3 HFPO-DA (IS)

331.8 &gt; 286.8 0.947 0.961 -0.014 671887 10.0 1479

1 Perfluoro(2-propoxypropanoic) acid

328.8 &gt; 284.8 0.961 0.988 -0.027 1.000 162499 2.17 56.3

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25017.d

Injection Date: 25-Jan-2018 11:14:03

Instrument ID: LC\_LCMS7

Lims ID: 280-105698-A-4-A

Lab Sample ID: 280-105698-4

Client ID: FAY-D-318BOONE-W1-1-012218

Operator ID: JBH

ALS Bottle#: 20

Worklist Smp#: 14

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

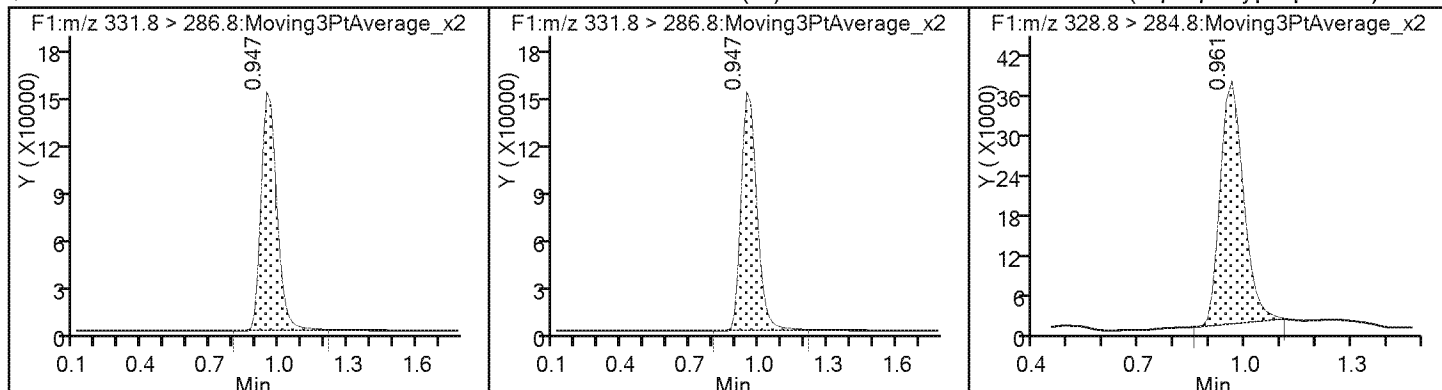

TestAmerica Denver  
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25017.d  
Lims ID: 280-105698-A-4-A  
Client ID: FAY-D-318BOONE-W1-1-012218  
Sample Type: Client  
Inject. Date: 25-Jan-2018 11:14:03 ALS Bottle#: 20 Worklist Smp#: 14  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: 280-105698-A-4-A  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:23 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyera

Date: 25-Jan-2018 15:28:34

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	9.19	91.86

FORM I  
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: TestAmerica Denver Job No.: 280-105698-1  
 SDG No.: \_\_\_\_\_  
 Client Sample ID: FAY-D-41BOONE-W1-1-012218 Lab Sample ID: 280-105698-5  
 Matrix: Water Lab File ID: hfpo718A25018.d  
 Analysis Method: 8321A Date Collected: 01/22/2018 13:49  
 Extraction Method: 3535 Date Extracted: 01/24/2018 15:20  
 Sample wt/vol: 252.7 (mL) Date Analyzed: 01/25/2018 11:17  
 Con. Extract Vol.: 5 (mL) Dilution Factor: 1  
 Injection Volume: 20 (uL) GC Column: Synergi Hydro ID: \_\_\_\_\_  
 % Moisture: \_\_\_\_\_ GPC Cleanup: (Y/N) N  
 Analysis Batch No.: 402806 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	<0.010		0.010	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	88		50-200



TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25018.d  
Lims ID: 280-105698-D-5-A  
Client ID: FAY-D-41BOONE-W1-1-012218  
Sample Type: Client  
Inject. Date: 25-Jan-2018 11:17:18 ALS Bottle#: 21 Worklist Smp#: 15  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: 280-105698-D-5-A  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:23 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyera

Date: 25-Jan-2018 15:28:36

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 &gt; 286.8 0.961 0.961 0.0 1.000 643864 8.80 1824

\* 2 13C3 HFPO-DA (IS)

331.8 &gt; 286.8 0.961 0.961 0.0 643864 10.0 1824

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25018.d

Injection Date: 25-Jan-2018 11:17:18

Instrument ID: LC\_LCMS7

Lims ID: 280-105698-D-5-A

Lab Sample ID: 280-105698-5

Client ID: FAY-D-41BOONE-W1-1-012218

Operator ID: JBH

ALS Bottle#: 21

Worklist Smp#: 15

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid (ND)

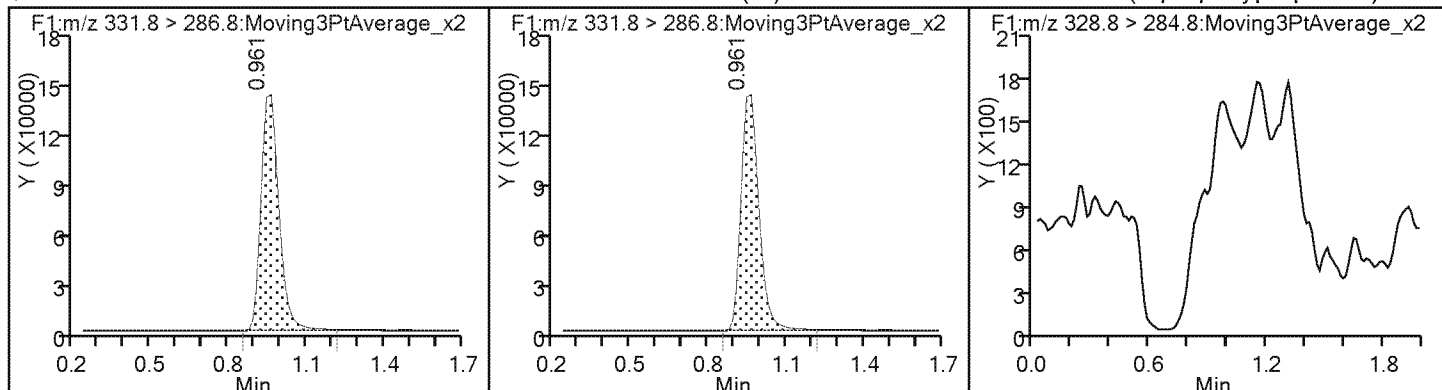

TestAmerica Denver  
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25018.d  
Lims ID: 280-105698-D-5-A  
Client ID: FAY-D-41BOONE-W1-1-012218  
Sample Type: Client  
Inject. Date: 25-Jan-2018 11:17:18 ALS Bottle#: 21 Worklist Smp#: 15  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: 280-105698-D-5-A  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:23 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyera

Date: 25-Jan-2018 15:28:36

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	8.80	88.03

FORM I  
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: <u>TestAmerica Denver</u>	Job No.: <u>280-105698-1</u>
SDG No.: _____	
Client Sample ID: <u>FAY-D-FB-012218</u>	Lab Sample ID: <u>280-105698-6</u>
Matrix: <u>Water</u>	Lab File ID: <u>hfpo718A25019.d</u>
Analysis Method: <u>8321A</u>	Date Collected: <u>01/22/2018 07:30</u>
Extraction Method: <u>3535</u>	Date Extracted: <u>01/24/2018 15:20</u>
Sample wt/vol: <u>260.4 (mL)</u>	Date Analyzed: <u>01/25/2018 11:20</u>
Con. Extract Vol.: <u>5 (mL)</u>	Dilution Factor: <u>1</u>
Injection Volume: <u>20 (uL)</u>	GC Column: <u>Synergi Hydro</u> ID: _____
% Moisture: _____	GPC Cleanup: (Y/N) <u>N</u>
Analysis Batch No.: <u>402806</u>	Units: <u>ug/L</u>

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	<0.010		0.010	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	92		50-200

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25019.d  
Lims ID: 280-105698-D-6-A  
Client ID: FAY-D-FB-012218  
Sample Type: Client  
Inject. Date: 25-Jan-2018 11:20:34 ALS Bottle#: 22 Worklist Smp#: 16  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: 280-105698-D-6-A  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:23 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyer

Date: 25-Jan-2018 15:28:39

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 &gt; 286.8 0.947 0.961 -0.014 1.000 674934 9.23 1466

\* 2 13C3 HFPO-DA (IS)

331.8 &gt; 286.8 0.947 0.961 -0.014 674934 10.0 1466

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25019.d

Injection Date: 25-Jan-2018 11:20:34

Instrument ID: LC\_LCMS7

Lims ID: 280-105698-D-6-A

Lab Sample ID: 280-105698-6

Client ID: FAY-D-FB-012218

Operator ID: JBH

ALS Bottle#: 22

Worklist Smp#: 16

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid (ND)

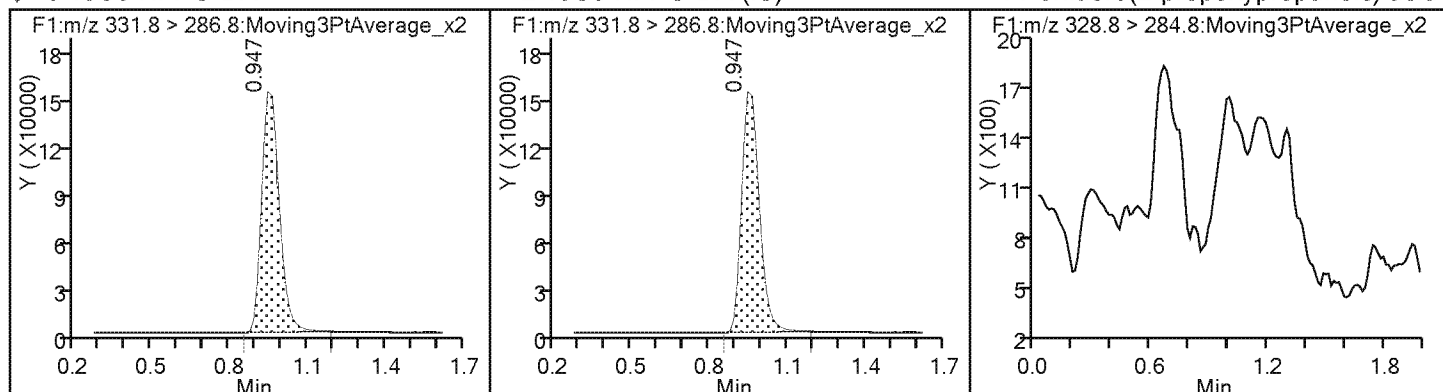

TestAmerica Denver  
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25019.d  
Lims ID: 280-105698-D-6-A  
Client ID: FAY-D-FB-012218  
Sample Type: Client  
Inject. Date: 25-Jan-2018 11:20:34 ALS Bottle#: 22 Worklist Smp#: 16  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: 280-105698-D-6-A  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:23 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyera

Date: 25-Jan-2018 15:28:39

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	9.23	92.27

Lab Name: <u>TestAmerica Denver</u>	Job No.: <u>280-105698-1</u>
SDG No.: _____	
Client Sample ID: <u>FAY-D-7145BUTLE-W1-1-0122</u> <u>18</u>	Lab Sample ID: <u>280-105698-7</u>
Matrix: <u>Water</u>	Lab File ID: <u>hfpo718A25020.d</u>
Analysis Method: <u>8321A</u>	Date Collected: <u>01/22/2018</u> <u>14:03</u>
Extraction Method: <u>3535</u>	Date Extracted: <u>01/24/2018</u> <u>15:20</u>
Sample wt/vol: <u>257.1(mL)</u>	Date Analyzed: <u>01/25/2018</u> <u>11:23</u>
Con. Extract Vol.: <u>5(mL)</u>	Dilution Factor: <u>1</u>
Injection Volume: <u>20(uL)</u>	GC Column: <u>Synergi Hydro</u> ID: _____
% Moisture: _____	GPC Cleanup: (Y/N) <u>N</u>
Analysis Batch No.: <u>402806</u>	Units: <u>ug/L</u>

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	0.080		0.010	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	84		50-200



TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25020.d  
Lims ID: 280-105698-C-7-A  
Client ID: FAY-D-7145BUTLE-W1-1-012218  
Sample Type: Client  
Inject. Date: 25-Jan-2018 11:23:49 ALS Bottle#: 23 Worklist Smp#: 17  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: 280-105698-C-7-A  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 02-Feb-2018 15:37:00 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK006

First Level Reviewer: meyera

Date: 25-Jan-2018 15:28:41

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 &gt; 286.8 0.947 0.961 -0.014 1.000 617530 8.44 1668

\* 2 13C3 HFPO-DA (IS)

331.8 &gt; 286.8 0.947 0.961 -0.014 617530 10.0 1668

1 Perfluoro(2-propoxypropanoic) acid

328.8 &gt; 284.8 0.947 0.988 -0.041 1.000 269901 4.10 74.0

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25020.d

Injection Date: 25-Jan-2018 11:23:49

Instrument ID: LC\_LCMS7

Lims ID: 280-105698-C-7-A

Lab Sample ID: 280-105698-7

Client ID: FAY-D-7145BUTLE-W1-1-012218

Operator ID: JBH

ALS Bottle#: 23

Worklist Smp#: 17

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

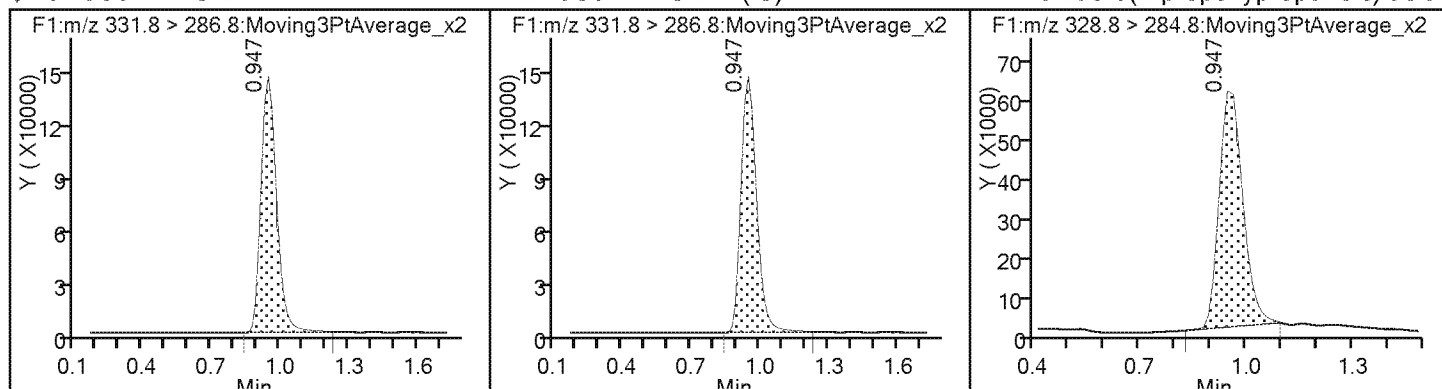

TestAmerica Denver  
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25020.d  
Lims ID: 280-105698-C-7-A  
Client ID: FAY-D-7145BUTLE-W1-1-012218  
Sample Type: Client  
Inject. Date: 25-Jan-2018 11:23:49 ALS Bottle#: 23 Worklist Smp#: 17  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: 280-105698-C-7-A  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 02-Feb-2018 15:37:00 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK006

First Level Reviewer: meyera

Date: 25-Jan-2018 15:28:41

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	8.44	84.43

## FORM I

Lab Name: <u>TestAmerica Denver</u>	Job No.: <u>280-105698-1</u>
SDG No.: _____	
Client Sample ID: <u>FAY-D-1515SCLLY-W1-1-0122</u> <u>18</u>	Lab Sample ID: <u>280-105698-8</u>
Matrix: <u>Water</u>	Lab File ID: <u>hfpo718A25021.d</u>
Analysis Method: <u>8321A</u>	Date Collected: <u>01/22/2018</u> <u>14:23</u>
Extraction Method: <u>3535</u>	Date Extracted: <u>01/24/2018</u> <u>15:20</u>
Sample wt/vol: <u>254.1(mL)</u>	Date Analyzed: <u>01/25/2018</u> <u>11:27</u>
Con. Extract Vol.: <u>5(mL)</u>	Dilution Factor: <u>1</u>
Injection Volume: <u>20(uL)</u>	GC Column: <u>Synergi Hydro</u> ID: _____
% Moisture: _____	GPC Cleanup: (Y/N) <u>N</u>
Analysis Batch No.: <u>402806</u>	Units: <u>ug/L</u>

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	<0.010		0.010	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	88		50-200

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25021.d  
Lims ID: 280-105698-B-8-A  
Client ID: FAY-D-1515SCLLY-W1-1-012218  
Sample Type: Client  
Inject. Date: 25-Jan-2018 11:27:05 ALS Bottle#: 24 Worklist Smp#: 18  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: 280-105698-B-8-A  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:23 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyer

Date: 25-Jan-2018 15:28:43

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 &gt; 286.8 0.961 0.961 0.0 1.000 641900 8.78 1425

\* 2 13C3 HFPO-DA (IS)

331.8 &gt; 286.8 0.961 0.961 0.0 641900 10.0 1425

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25021.d

Injection Date: 25-Jan-2018 11:27:05

Instrument ID: LC\_LCMS7

Lims ID: 280-105698-B-8-A

Lab Sample ID: 280-105698-8

Client ID: FAY-D-1515SCLLY-W1-1-012218

Operator ID: JBH

ALS Bottle#: 24

Worklist Smp#: 18

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid (ND)

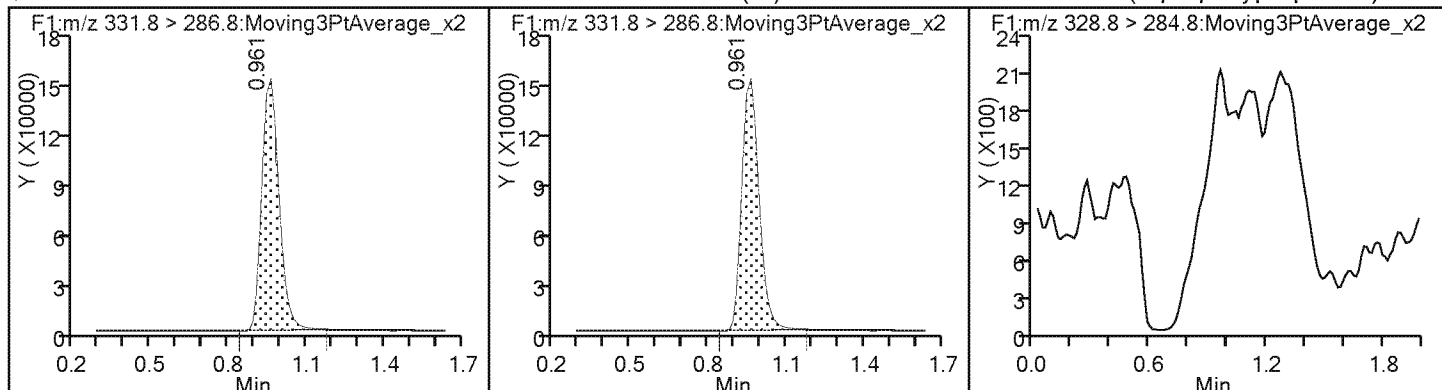

TestAmerica Denver  
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25021.d  
Lims ID: 280-105698-B-8-A  
Client ID: FAY-D-1515SCLLY-W1-1-012218  
Sample Type: Client  
Inject. Date: 25-Jan-2018 11:27:05 ALS Bottle#: 24 Worklist Smp#: 18  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: 280-105698-B-8-A  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:23 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyera

Date: 25-Jan-2018 15:28:43

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	8.78	87.76





TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25023.d  
Lims ID: 280-105698-A-9-A  
Client ID: FAY-D-7396SALIE-W1-1-012218  
Sample Type: Client  
Inject. Date: 25-Jan-2018 11:33:36 ALS Bottle#: 25 Worklist Smp#: 20  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: 280-105698-A-9-A  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:31 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyera

Date: 25-Jan-2018 15:28:48

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 &gt; 286.8 1.002 0.961 0.041 1.000 630863 8.62 1381

\* 2 13C3 HFPO-DA (IS)

331.8 &gt; 286.8 1.002 0.961 0.041 630863 10.0 1381

1 Perfluoro(2-propoxypropanoic) acid

328.8 &gt; 284.8 1.015 0.988 0.027 1.000 28439 0.2295 7.4

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25023.d

Injection Date: 25-Jan-2018 11:33:36

Instrument ID: LC\_LCMS7

Lims ID: 280-105698-A-9-A

Lab Sample ID: 280-105698-9

Client ID: FAY-D-7396SALIE-W1-1-012218

Operator ID: JBH

ALS Bottle#: 25

Worklist Smp#: 20

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

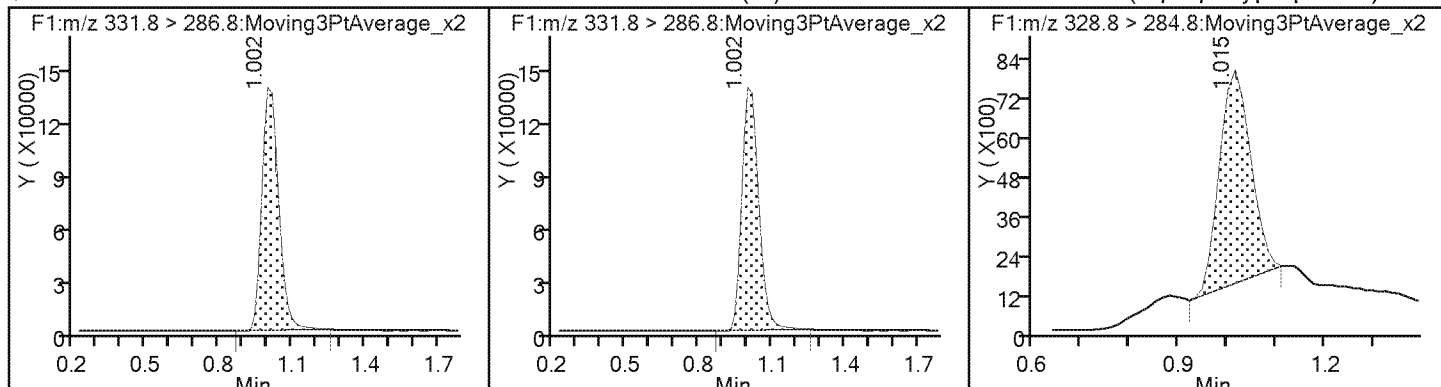

TestAmerica Denver  
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25023.d  
Lims ID: 280-105698-A-9-A  
Client ID: FAY-D-7396SALIE-W1-1-012218  
Sample Type: Client  
Inject. Date: 25-Jan-2018 11:33:36 ALS Bottle#: 25 Worklist Smp#: 20  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: 280-105698-A-9-A  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:31 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyera

Date: 25-Jan-2018 15:28:48

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	8.62	86.25

Lab Name: <u>TestAmerica Denver</u>	Job No.: <u>280-105698-1</u>
SDG No.: _____	
Client Sample ID: <u>FAY-D-7012NC87H-W1-1-0122</u> <u>18</u>	Lab Sample ID: <u>280-105698-10</u>
Matrix: <u>Water</u>	Lab File ID: <u>hfpo718A25024.d</u>
Analysis Method: <u>8321A</u>	Date Collected: <u>01/22/2018</u> <u>16:21</u>
Extraction Method: <u>3535</u>	Date Extracted: <u>01/24/2018</u> <u>15:20</u>
Sample wt/vol: <u>247.7(mL)</u>	Date Analyzed: <u>01/25/2018</u> <u>11:36</u>
Con. Extract Vol.: <u>5(mL)</u>	Dilution Factor: <u>1</u>
Injection Volume: <u>20(uL)</u>	GC Column: <u>Synergi Hydro</u> ID: _____
% Moisture: _____	GPC Cleanup: (Y/N) <u>N</u>
Analysis Batch No.: <u>402806</u>	Units: <u>ug/L</u>

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	0.027		0.010	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	86		50-200

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25024.d  
Lims ID: 280-105698-A-10-A  
Client ID: FAY-D-7012NC87-W1-1-012218  
Sample Type: Client  
Inject. Date: 25-Jan-2018 11:36:52 ALS Bottle#: 26 Worklist Smp#: 21  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: 280-105698-A-10-A  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:31 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyer

Date: 25-Jan-2018 15:28:51

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 &gt; 286.8 0.975 0.961 0.014 1.000 629975 8.61 1450

\* 2 13C3 HFPO-DA (IS)

331.8 &gt; 286.8 0.975 0.961 0.014 629975 10.0 1450

1 Perfluoro(2-propoxypropanoic) acid

328.8 &gt; 284.8 0.975 0.988 -0.013 1.000 98892 1.34 42.8

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25024.d

Injection Date: 25-Jan-2018 11:36:52

Instrument ID: LC\_LCMS7

Lims ID: 280-105698-A-10-A

Lab Sample ID: 280-105698-10

Client ID: FAY-D-7012NC87-W1-1-012218

Operator ID: JBH

ALS Bottle#: 26

Worklist Smp#: 21

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

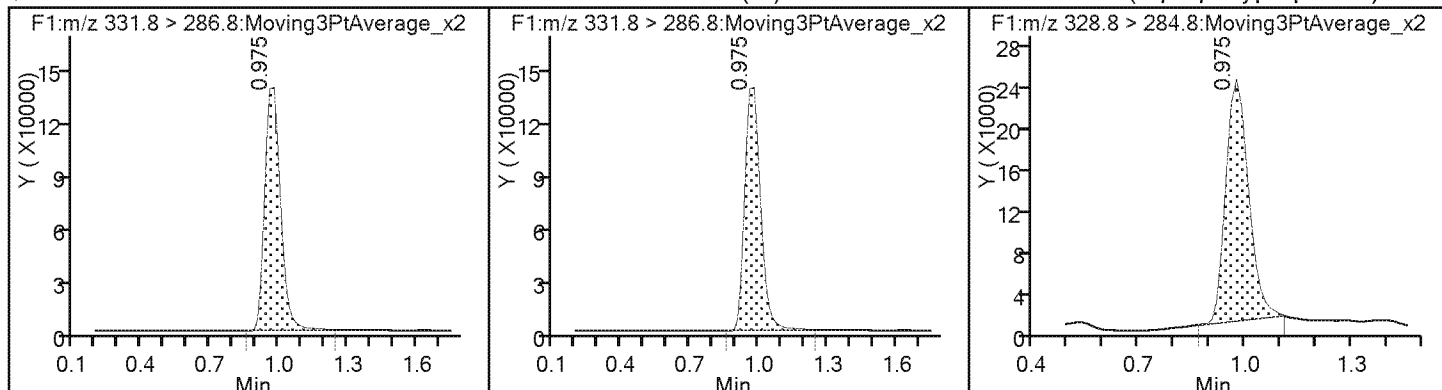

TestAmerica Denver  
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25024.d  
Lims ID: 280-105698-A-10-A  
Client ID: FAY-D-7012NC87-W1-1-012218  
Sample Type: Client  
Inject. Date: 25-Jan-2018 11:36:52 ALS Bottle#: 26 Worklist Smp#: 21  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: 280-105698-A-10-A  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:31 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyera

Date: 25-Jan-2018 15:28:51

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	8.61	86.13

FORM VI  
LCMS BY EXTERNAL STANDARD - INITIAL CALIBRATION DATA  
RETENTION TIME SUMMARY

Lab Name: TestAmerica Denver Job No.: 280-105698-1 Analy Batch No.: 387775  
 SDG No.: \_\_\_\_\_  
 Instrument ID: LC\_LCMS7 GC Column: Synergi Hyd ID: \_\_\_\_\_ Heated Purge: (Y/N) N  
 Calibration Start Date: 09/14/2017 14:40 Calibration End Date: 09/14/2017 15:01 Calibration ID: 30321

Calibration Files:

LEVEL:	LAB SAMPLE ID:	LAB FILE ID:
Level 1	STD001 280-387775/3	hfpo717I14052.d
Level 2	STD002 280-387775/4	hfpo717I14053.d
Level 3	STD003 280-387775/5	hfpo717I14054.d
Level 4	STD004 280-387775/6	hfpo717I14055.d
Level 5	STD005 280-387775/7	hfpo717I14056.d
Level 6	STD006 280-387775/8	hfpo717I14057.d
Level 7	STD007 280-387775/9	hfpo717I14058.d
Level 8	STD008 280-387775/10	hfpo717I14059.d

ANALYTE	LVL 1	LVL 2	LVL 3	LVL 4	LVL 5	LVL 6	LVL 7	LVL 8			RT WINDOW	AVG RT
Perfluoro(2-propoxypropanoic) acid	1.002	0.988	0.988	0.988	0.975	0.975	0.988	0.988			0.486 - 1.486	0.987
13C3 HFPO-DA	0.988	0.975	0.975	0.988	0.975	0.975	0.988	0.988			0.481 - 1.481	0.982



FORM VI  
LCMS BY EXTERNAL STANDARD - INITIAL CALIBRATION DATA  
CURVE EVALUATION

Lab Name: TestAmerica Denver Job No.: 280-105698-1 Analy Batch No.: 387775  
 SDG No.: \_\_\_\_\_  
 Instrument ID: LC\_LCMS7 GC Column: Synergi Hyd ID: \_\_\_\_\_ Heated Purge: (Y/N) N  
 Calibration Start Date: 09/14/2017 14:40 Calibration End Date: 09/14/2017 15:01 Calibration ID: 30321

Calibration Files:

LEVEL:	LAB SAMPLE ID:	LAB FILE ID:
Level 1	STD001 280-387775/3	hfpo717I14052.d
Level 2	STD002 280-387775/4	hfpo717I14053.d
Level 3	STD003 280-387775/5	hfpo717I14054.d
Level 4	STD004 280-387775/6	hfpo717I14055.d
Level 5	STD005 280-387775/7	hfpo717I14056.d
Level 6	STD006 280-387775/8	hfpo717I14057.d
Level 7	STD007 280-387775/9	hfpo717I14058.d
Level 8	STD008 280-387775/10	hfpo717I14059.d

ANALYTE	CF				CURVE TYPE	COEFFICIENT			#	MIN CF	%RSD	#	MAX %RSD	R^2 OR COD	#	MIN R^2 OR COD
	LVL 1 LVL 5	LVL 2 LVL 6	LVL 3 LVL 7	LVL 4 LVL 8		B	M1	M2								
13C3 HFPO-DA	206978 182499	200375 172499	208177 188495	195084 187811	Ave		192739.525				6.4		30.0			

Note: The M1 coefficient is the same as Ave CF for an Ave curve type.

FORM VI  
LCMS BY INTERNAL STANDARD - INITIAL CALIBRATION DATA  
CURVE EVALUATION

Lab Name: TestAmerica Denver Job No.: 280-105698-1 Analy Batch No.: 387775  
 SDG No.: \_\_\_\_\_  
 Instrument ID: LC\_LCMS7 GC Column: Synergi Hyd ID: \_\_\_\_\_ Heated Purge: (Y/N) N  
 Calibration Start Date: 09/14/2017 14:40 Calibration End Date: 09/14/2017 15:01 Calibration ID: 30321

ANALYTE	RRF					CURVE TYPE	COEFFICIENT			#	MIN RRF	%RSD	#	MAX %RSD	R^2 OR COD	#	MIN R^2 OR COD
	LVL 1 LVL 6	LVL 2 LVL 7	LVL 3 LVL 8	LVL 4	LVL 5		B	M1	M2								
Perfluoro(2-propoxypropanoic) acid	1.6312 0.9640	1.1780 0.9353	0.9745 0.8831	0.9868	1.0688	Lin1	0.1732	0.9076							0.9980		0.9900

Note: The M1 coefficient is the same as Ave RRF for an Ave curve type.

FORM VI  
LCMS BY EXTERNAL STANDARD - INITIAL CALIBRATION DATA  
RESPONSE AND CONCENTRATION

Lab Name: TestAmerica Denver Job No.: 280-105698-1 Analy Batch No.: 387775  
SDG No.: \_\_\_\_\_  
Instrument ID: LC\_LCMS7 GC Column: Synergi Hyd ID: \_\_\_\_\_ Heated Purge: (Y/N) N  
Calibration Start Date: 09/14/2017 14:40 Calibration End Date: 09/14/2017 15:01 Calibration ID: 30321

Calibration Files:

LEVEL:	LAB SAMPLE ID:	LAB FILE ID:
Level 1	STD001 280-387775/3	hfpo717I14052.d
Level 2	STD002 280-387775/4	hfpo717I14053.d
Level 3	STD003 280-387775/5	hfpo717I14054.d
Level 4	STD004 280-387775/6	hfpo717I14055.d
Level 5	STD005 280-387775/7	hfpo717I14056.d
Level 6	STD006 280-387775/8	hfpo717I14057.d
Level 7	STD007 280-387775/9	hfpo717I14058.d
Level 8	STD008 280-387775/10	hfpo717I14059.d

ANALYTE	CURVE TYPE	RESPONSE					CONCENTRATION (UG/L)				
		LVL 1 LVL 6	LVL 2 LVL 7	LVL 3 LVL 8	LVL 4	LVL 5	LVL 1 LVL 6	LVL 2 LVL 7	LVL 3 LVL 8	LVL 4	LVL 5
13C3 HFPO-DA	Ave	2069777	2003748	2081766	1950837	1824991	10.0	10.0	10.0	10.0	10.0
		1724989	1884947	1878107			10.0	10.0	10.0		

Curve Type Legend:

Ave = Average

FORM VI  
LCMS BY INTERNAL STANDARD - INITIAL CALIBRATION DATA  
RESPONSE AND CONCENTRATION

Lab Name: TestAmerica Denver Job No.: 280-105698-1 Analy Batch No.: 387775

SDG No.: \_\_\_\_\_

Instrument ID: LC\_LCMS7 GC Column: Synergi Hyd ID: \_\_\_\_\_ Heated Purge: (Y/N) N

Calibration Start Date: 09/14/2017 14:40 Calibration End Date: 09/14/2017 15:01 Calibration ID: 30321

Calibration Files:

LEVEL:	LAB SAMPLE ID:	LAB FILE ID:
Level 1	STD001 280-387775/3	hfpo717I14052.d
Level 2	STD002 280-387775/4	hfpo717I14053.d
Level 3	STD003 280-387775/5	hfpo717I14054.d
Level 4	STD004 280-387775/6	hfpo717I14055.d
Level 5	STD005 280-387775/7	hfpo717I14056.d
Level 6	STD006 280-387775/8	hfpo717I14057.d
Level 7	STD007 280-387775/9	hfpo717I14058.d
Level 8	STD008 280-387775/10	hfpo717I14059.d

ANALYTE	IS REF	CURVE TYPE	RESPONSE					CONCENTRATION (UG/L)				
			LVL 1 LVL 6	LVL 2 LVL 7	LVL 3 LVL 8	LVL 4	LVL 5	LVL 1 LVL 6	LVL 2 LVL 7	LVL 3 LVL 8	LVL 4	LVL 5
Perfluoro(2-propoxypropanoic) acid	13CP ODA	Lin1	84406	118017	202876	385009	975278	0.250	0.500	1.00	2.00	5.00
			1662919	4407541	8293101			10.0	25.0	50.0		

Curve Type Legend:

Lin1 = Linear 1/conc ISTD

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\hfpo717114052.d  
Lims ID: std001  
Client ID:  
Sample Type: IC Calib Level: 1  
Inject. Date: 14-Sep-2017 14:40:03 ALS Bottle#: 2 Worklist Smp#: 3  
Injection Vol: 10.0 ul Dil. Factor: 1.0000  
Sample Info: L1  
Misc. Info.: HFPO17114  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Sublist: chrom-HFPO\*sub1  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 15-Sep-2017 07:29:39 Calib Date: 14-Sep-2017 15:01:22  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\hfpo717114059.d

Column 1 : Det: F1:MRM  
Process Host: XAWRK034

First Level Reviewer: meyera Date: 15-Sep-2017 07:28:15

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 > 286.8 0.988 0.981 0.007 1.000 2069777 10.7 429

\* 2 13C3 HFPO-DA (IS)

331.8 > 286.8 0.988 0.981 0.007 2069777 10.0 429

1 Perfluoro(2-propoxypropanoic) acid

328.8 > 284.8 1.002 0.986 0.016 1.000 84406 0.2585 49.7

**Reagents:**

HFPO\_CAL-1\_00030

Amount Added: 1.00

Units: mL

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\hfpo717114052.d

Injection Date: 14-Sep-2017 14:40:03

Instrument ID: LC\_LCMS7

Lims ID: std001

Client ID:

Operator ID: JBH

ALS Bottle#:

2

Worklist Smp#:

3

Injection Vol: 10.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

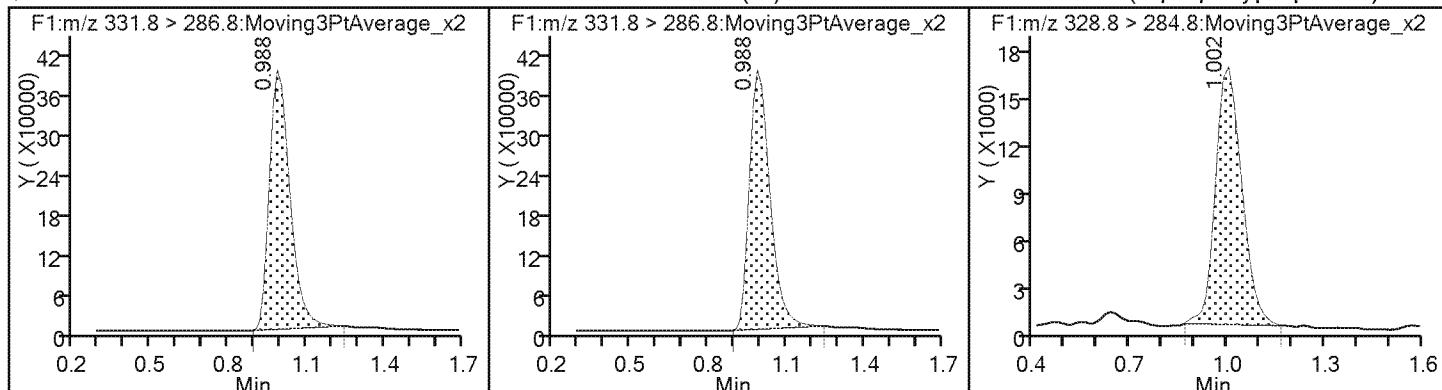

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\hfpo717114053.d  
Lims ID: std002  
Client ID:  
Sample Type: IC Calib Level: 2  
Inject. Date: 14-Sep-2017 14:43:06 ALS Bottle#: 3 Worklist Smp#: 4  
Injection Vol: 10.0 ul Dil. Factor: 1.0000  
Sample Info: L2  
Misc. Info.: HFPO17114  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Sublist: chrom-HFPO\*sub1  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 15-Sep-2017 07:29:39 Calib Date: 14-Sep-2017 15:01:22  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\hfpo717114059.d

Column 1 : Det: F1:MRM  
Process Host: XAWRK034

First Level Reviewer: meyera Date: 15-Sep-2017 07:28:18

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\* 2 13C3 HFPO-DA (IS)

331.8 > 286.8 0.975 0.981 -0.006 2003748 10.0 386

\$ 3 13C3 HFPO-DA

331.8 > 286.8 0.975 0.981 -0.006 1.000 2003748 10.4 386

1 Perfluoro(2-propoxypropanoic) acid

328.8 > 284.8 0.988 0.986 0.002 1.000 118017 0.4581 56.6

**Reagents:**

HFPO\_CAL-2\_00031

Amount Added: 1.00

Units: mL

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\hfpo717114053.d

Injection Date: 14-Sep-2017 14:43:06

Instrument ID: LC\_LCMS7

Lims ID: std002

Client ID:

Operator ID: JBH

ALS Bottle#: 3

Worklist Smp#: 4

Injection Vol: 10.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\* 2 13C3 HFPO-DA (IS)

\$ 3 13C3 HFPO-DA

1 Perfluoro(2-propoxypropanoic) acid

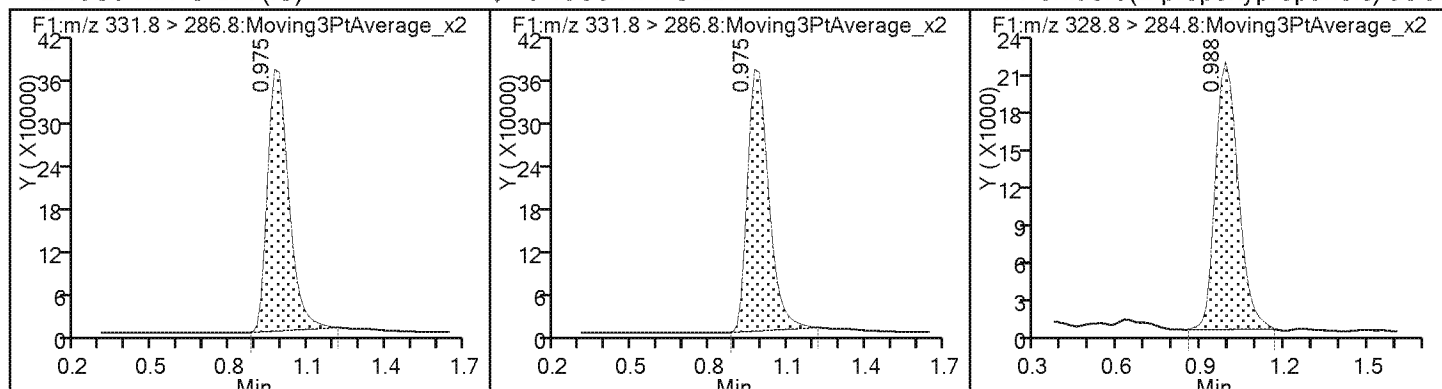



TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\hfpo717114054.d  
 Lims ID: std003  
 Client ID:  
 Sample Type: IC Calib Level: 3  
 Inject. Date: 14-Sep-2017 14:46:08 ALS Bottle#: 4 Worklist Smp#: 5  
 Injection Vol: 10.0 ul Dil. Factor: 1.0000  
 Sample Info: L3  
 Misc. Info.: HFPO17114  
 Operator ID: JBH Instrument ID: LC\_LCMS7  
 Sublist: chrom-HFPO\*sub1  
 Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\HFPO.m  
 Limit Group: LC - 8321A\_HFPO\_Du  
 Last Update: 15-Sep-2017 07:29:40 Calib Date: 14-Sep-2017 15:01:22  
 Integrator: Picker  
 Quant Method: Internal/External Standard Quant By: Initial Calibration  
 Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\hfpo717114059.d

Column 1 : Det: F1:MRM  
 Process Host: XAWRK034

First Level Reviewer: meyer Date: 15-Sep-2017 07:28:20

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 > 286.8 0.975 0.981 -0.006 1.000 2081766 10.8 403

\* 2 13C3 HFPO-DA (IS)

331.8 > 286.8 0.975 0.981 -0.006 2081766 10.0 403

1 Perfluoro(2-propoxypropanoic) acid

328.8 > 284.8 0.988 0.986 0.002 1.000 202876 0.8830 108

**Reagents:**

HFPO\_CAL-3\_00030 Amount Added: 1.00 Units: mL

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\hfpo717114054.d

Injection Date: 14-Sep-2017 14:46:08

Instrument ID: LC\_LCMS7

Lims ID: std003

Client ID:

Operator ID: JBH

ALS Bottle#:

4

Worklist Smp#:

5

Injection Vol: 10.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

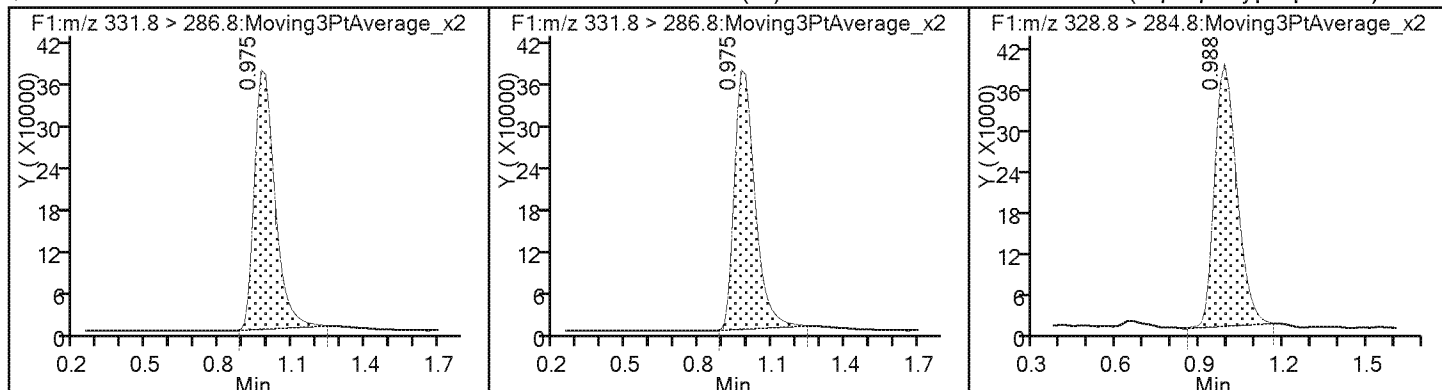

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\hfpo717114055.d  
Lims ID: std004  
Client ID:  
Sample Type: IC Calib Level: 4  
Inject. Date: 14-Sep-2017 14:49:11 ALS Bottle#: 5 Worklist Smp#: 6  
Injection Vol: 10.0 ul Dil. Factor: 1.0000  
Sample Info: L4  
Misc. Info.: HFPO17114  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Sublist: chrom-HFPO\*sub1  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 15-Sep-2017 07:29:40 Calib Date: 14-Sep-2017 15:01:22  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\hfpo717114059.d

Column 1 : Det: F1:MRM  
Process Host: XAWRK034

First Level Reviewer: meyera Date: 15-Sep-2017 07:28:22

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\* 2 13C3 HFPO-DA (IS)

331.8 > 286.8 0.988 0.981 0.007 1950837 10.0 384

\$ 3 13C3 HFPO-DA

331.8 > 286.8 0.988 0.981 0.007 1.000 1950837 10.1 384

1 Perfluoro(2-propoxypropanoic) acid

328.8 > 284.8 0.988 0.986 0.002 1.000 385009 1.98 162

**Reagents:**

HFPO\_CAL-4\_00030 Amount Added: 1.00 Units: mL

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\hfp0717114055.d

Injection Date: 14-Sep-2017 14:49:11

Instrument ID: LC\_LCMS7

Lims ID: std004

Client ID:

Operator ID: JBH

ALS Bottle#: 5

Worklist Smp#: 6

Injection Vol: 10.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\* 2 13C3 HFPO-DA (IS)

\$ 3 13C3 HFPO-DA

1 Perfluoro(2-propoxypropanoic) acid

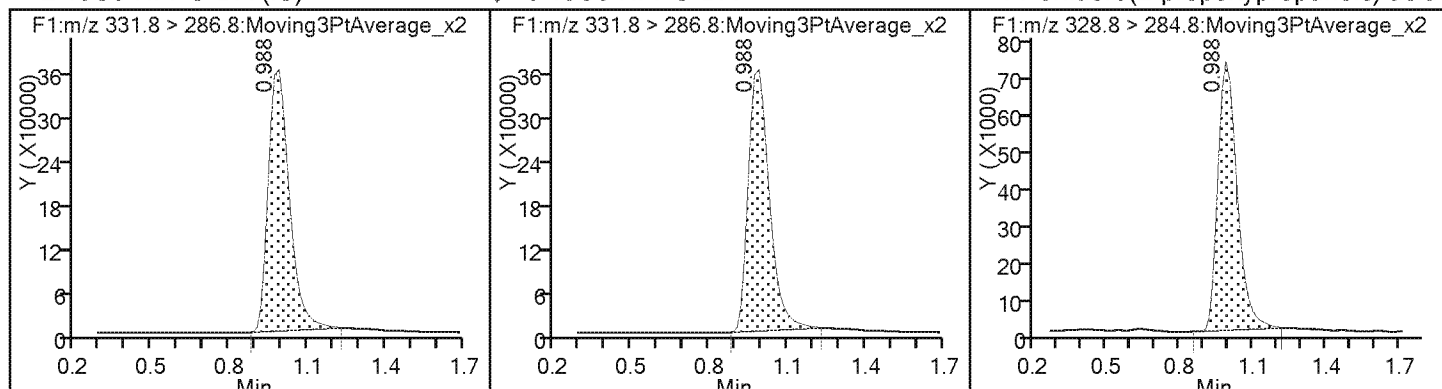

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\hfpo717114056.d  
Lims ID: std005  
Client ID:  
Sample Type: IC Calib Level: 5  
Inject. Date: 14-Sep-2017 14:52:13 ALS Bottle#: 6 Worklist Smp#: 7  
Injection Vol: 10.0 ul Dil. Factor: 1.0000  
Sample Info: L5  
Misc. Info.: HFPO17114  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Sublist: chrom-HFPO\*sub1  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 15-Sep-2017 07:29:41 Calib Date: 14-Sep-2017 15:01:22  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\hfpo717114059.d

Column 1 : Det: F1:MRM  
Process Host: XAWRK034

First Level Reviewer: meyera Date: 15-Sep-2017 07:28:25

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 > 286.8 0.975 0.981 -0.006 1.000 1824991 9.47 371

\* 2 13C3 HFPO-DA (IS)

331.8 > 286.8 0.975 0.981 -0.006 1824991 10.0 371

1 Perfluoro(2-propoxypropanoic) acid

328.8 > 284.8 0.975 0.986 -0.011 1.000 975278 5.70 268

**Reagents:**

HFPO\_CAL-5\_00067 Amount Added: 1.00 Units: mL

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\hfpo717114056.d

Injection Date: 14-Sep-2017 14:52:13

Instrument ID: LC\_LCMS7

Lims ID: std005

Client ID:

Operator ID: JBH

ALS Bottle#: 6

Worklist Smp#: 7

Injection Vol: 10.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

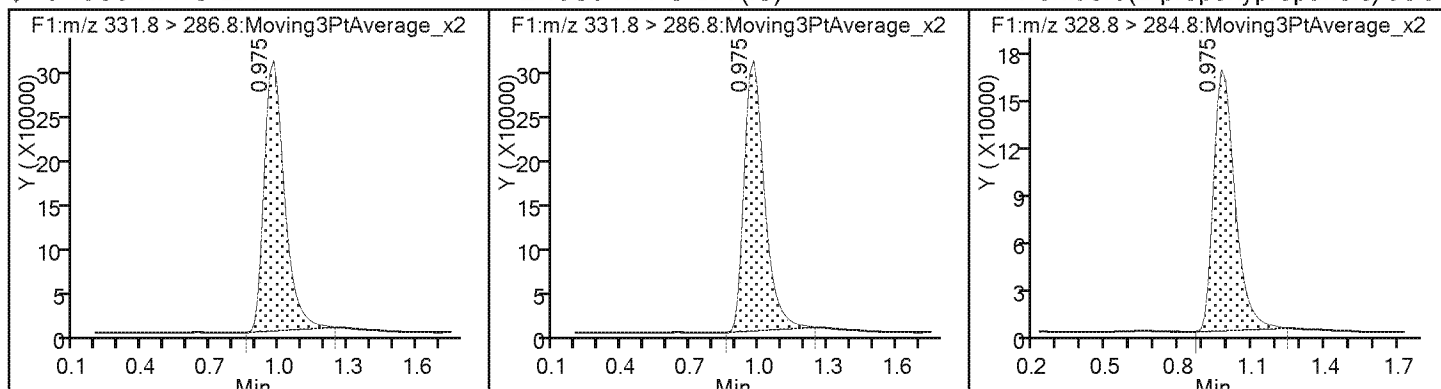

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\hfpo717114057.d  
Lims ID: std006  
Client ID:  
Sample Type: IC Calib Level: 6  
Inject. Date: 14-Sep-2017 14:55:16 ALS Bottle#: 7 Worklist Smp#: 8  
Injection Vol: 10.0 ul Dil. Factor: 1.0000  
Sample Info: L6  
Misc. Info.: HFPO17114  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Sublist: chrom-HFPO\*sub1  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 15-Sep-2017 07:29:41 Calib Date: 14-Sep-2017 15:01:22  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\hfpo717114059.d

Column 1 : Det: F1:MRM  
Process Host: XAWRK034

First Level Reviewer: meyera Date: 15-Sep-2017 07:28:27

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\* 2 13C3 HFPO-DA (IS)

331.8 > 286.8 0.975 0.981 -0.006 1724989 10.0 287

\$ 3 13C3 HFPO-DA

331.8 > 286.8 0.975 0.981 -0.006 1.000 1724989 8.95 287

1 Perfluoro(2-propoxypropanoic) acid

328.8 > 284.8 0.975 0.986 -0.011 1.000 1662919 10.4 248

**Reagents:**

HFPO\_CAL-6\_00067 Amount Added: 1.00 Units: mL

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\hfpo717114057.d

Injection Date: 14-Sep-2017 14:55:16

Instrument ID: LC\_LCMS7

Lims ID: std006

Client ID:

Operator ID: JBH

ALS Bottle#: 7

Worklist Smp#: 8

Injection Vol: 10.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\* 2 13C3 HFPO-DA (IS)

\$ 3 13C3 HFPO-DA

1 Perfluoro(2-propoxypropanoic) acid

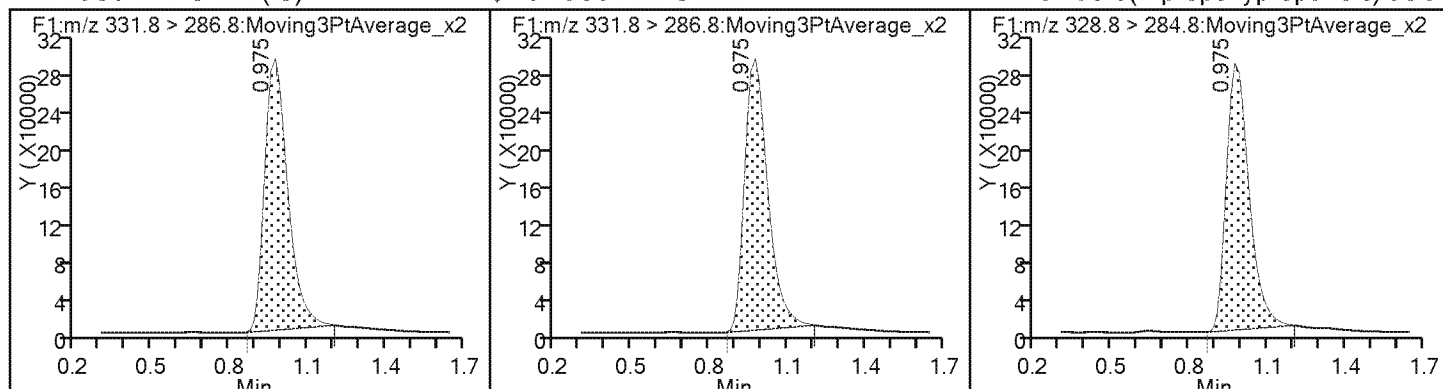



TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\hfpo717114058.d  
Lims ID: std007  
Client ID:  
Sample Type: IC Calib Level: 7  
Inject. Date: 14-Sep-2017 14:58:19 ALS Bottle#: 8 Worklist Smp#: 9  
Injection Vol: 10.0 ul Dil. Factor: 1.0000  
Sample Info: L7  
Misc. Info.: HFPO17114  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Sublist: chrom-HFPO\*sub1  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 15-Sep-2017 07:29:42 Calib Date: 14-Sep-2017 15:01:22  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\hfpo717114059.d

Column 1 : Det: F1:MRM  
Process Host: XAWRK034

First Level Reviewer: meyer Date: 15-Sep-2017 07:28:29

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 > 286.8 0.988 0.981 0.007 1.000 1884947 9.78 361

\* 2 13C3 HFPO-DA (IS)

331.8 > 286.8 0.988 0.981 0.007 1884947 10.0 361

1 Perfluoro(2-propoxypropanoic) acid

328.8 > 284.8 0.988 0.986 0.002 1.000 4407541 25.6 379

**Reagents:**

HFPO\_CAL-7\_00030

Amount Added: 1.00

Units: mL

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\hfpo717114058.d

Injection Date: 14-Sep-2017 14:58:19

Instrument ID: LC\_LCMS7

Lims ID: std007

Client ID:

Operator ID: JBH

ALS Bottle#: 8

Worklist Smp#: 9

Injection Vol: 10.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

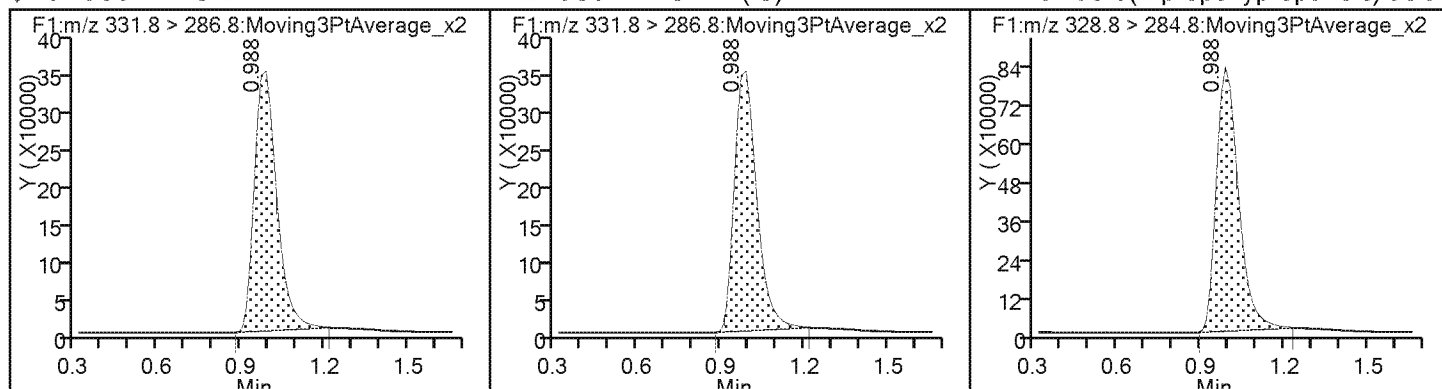

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\hfpo717114059.d  
Lims ID: std008  
Client ID:  
Sample Type: IC Calib Level: 8  
Inject. Date: 14-Sep-2017 15:01:22 ALS Bottle#: 9 Worklist Smp#: 10  
Injection Vol: 10.0 ul Dil. Factor: 1.0000  
Sample Info: L8  
Misc. Info.: HFPO17114  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Sublist: chrom-HFPO\*sub1  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 15-Sep-2017 07:29:43 Calib Date: 14-Sep-2017 15:01:22  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\hfpo717114059.d

Column 1 : Det: F1:MRM  
Process Host: XAWRK034

First Level Reviewer: meyera Date: 15-Sep-2017 07:28:32

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\* 2 13C3 HFPO-DA (IS)

331.8 > 286.8 0.988 0.981 0.007 1878107 10.0 379

\$ 3 13C3 HFPO-DA

331.8 > 286.8 0.988 0.981 0.007 1.000 1878107 9.74 379

1 Perfluoro(2-propoxypropanoic) acid

328.8 > 284.8 0.988 0.986 0.002 1.000 8293101 48.5 359

**Reagents:**

HFPO\_CAL-8\_00030

Amount Added: 1.00

Units: mL

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\hfpo717114059.d

Injection Date: 14-Sep-2017 15:01:22

Instrument ID: LC\_LCMS7

Lims ID: std008

Client ID:

Operator ID: JBH

ALS Bottle#:

9

Worklist Smp#:

10

Injection Vol: 10.0 ul


Dil. Factor: 1.0000

Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\* 2 <sup>13</sup>C3 HFPO-DA (IS)\$ 3 <sup>13</sup>C3 HFPO-DA

1 Perfluoro(2-propoxypropanoic) acid

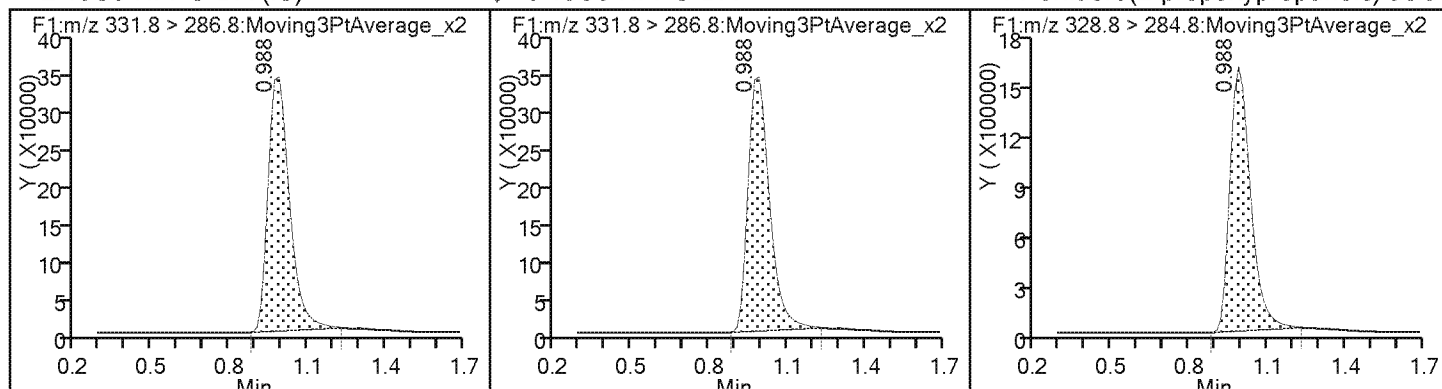

FORM VI  
LCMS BY EXTERNAL STANDARD - INITIAL CALIBRATION DATA  
RETENTION TIME SUMMARY

Lab Name: TestAmerica Denver Job No.: 280-105698-1 Analy Batch No.: 390728  
SDG No.: \_\_\_\_\_  
Instrument ID: LC\_LCMS7 GC Column: Synergi Hyd ID: \_\_\_\_\_ Heated Purge: (Y/N) N  
Calibration Start Date: 10/10/2017 09:35 Calibration End Date: 10/10/2017 09:58 Calibration ID: 30558

Calibration Files:

LEVEL:	LAB SAMPLE ID:	LAB FILE ID:
Level 1	STD001 280-390728/3	hfpo717J10026.d
Level 2	STD002 280-390728/4	hfpo717J10027.d
Level 3	STD003 280-390728/5	hfpo717J10028.d
Level 4	STD004 280-390728/6	hfpo717J10029.d
Level 5	STD005 280-390728/7	hfpo717J10030.d
Level 6	STD006 280-390728/8	hfpo717J10031.d
Level 7	STD007 280-390728/9	hfpo717J10032.d
Level 8	STD008 280-390728/10	hfpo717J10033.d

ANALYTE	LVL 1	LVL 2	LVL 3	LVL 4	LVL 5	LVL 6	LVL 7	LVL 8			RT WINDOW	AVG RT
HFPO-DA	0.893	0.880	0.880	0.880	0.893	0.880	0.880	0.893			0.385 - 1.385	0.885
13C3 HFPO-DA	0.880	0.880	0.880	0.880	0.880	0.880	0.880	0.880			0.380 - 1.380	0.880

FORM VI  
LCMS BY EXTERNAL STANDARD - INITIAL CALIBRATION DATA  
CURVE EVALUATION

Lab Name: TestAmerica Denver Job No.: 280-105698-1 Analy Batch No.: 390728  
 SDG No.: \_\_\_\_\_  
 Instrument ID: LC\_LCMS7 GC Column: Synergi Hyd ID: \_\_\_\_\_ Heated Purge: (Y/N) N  
 Calibration Start Date: 10/10/2017 09:35 Calibration End Date: 10/10/2017 09:58 Calibration ID: 30558

Calibration Files:

LEVEL:	LAB SAMPLE ID:	LAB FILE ID:
Level 1	STD001 280-390728/3	hfpo717J10026.d
Level 2	STD002 280-390728/4	hfpo717J10027.d
Level 3	STD003 280-390728/5	hfpo717J10028.d
Level 4	STD004 280-390728/6	hfpo717J10029.d
Level 5	STD005 280-390728/7	hfpo717J10030.d
Level 6	STD006 280-390728/8	hfpo717J10031.d
Level 7	STD007 280-390728/9	hfpo717J10032.d
Level 8	STD008 280-390728/10	hfpo717J10033.d

ANALYTE	CF				CURVE TYPE	COEFFICIENT			#	MIN CF	%RSD	#	MAX %RSD	R^2 OR COD	#	MIN R^2 OR COD
	LVL 1 LVL 5	LVL 2 LVL 6	LVL 3 LVL 7	LVL 4 LVL 8		B	M1	M2								
13C3 HFPO-DA	73075 74460	74523 73194	75043 72919	71803 70142	Ave		73144.6750				2.2		30.0			

Note: The M1 coefficient is the same as Ave CF for an Ave curve type.

FORM VI  
LCMS BY INTERNAL STANDARD - INITIAL CALIBRATION DATA  
CURVE EVALUATION

Lab Name: TestAmerica Denver Job No.: 280-105698-1 Analy Batch No.: 390728  
 SDG No.: \_\_\_\_\_  
 Instrument ID: LC\_LCMS7 GC Column: Synergi Hyd ID: \_\_\_\_\_ Heated Purge: (Y/N) N  
 Calibration Start Date: 10/10/2017 09:35 Calibration End Date: 10/10/2017 09:58 Calibration ID: 30558

ANALYTE	RRF					CURVE TYPE	COEFFICIENT			#	MIN RRF	%RSD	#	MAX %RSD	R^2 OR COD	#	MIN R^2 OR COD
	LVL 1 LVL 6	LVL 2 LVL 7	LVL 3 LVL 8	LVL 4	LVL 5		B	M1	M2								
HFPO-DA	1.6980 1.0102	1.7128 0.9824	1.1896 1.0419	1.1637	1.0154	Lin1	0.2185	1.0121							0.9980		0.9900

Note: The M1 coefficient is the same as Ave RRF for an Ave curve type.

FORM VI  
LCMS BY EXTERNAL STANDARD - INITIAL CALIBRATION DATA  
RESPONSE AND CONCENTRATION

Lab Name: TestAmerica Denver Job No.: 280-105698-1 Analy Batch No.: 390728  
 SDG No.: \_\_\_\_\_  
 Instrument ID: LC\_LCMS7 GC Column: Synergi Hyd ID: \_\_\_\_\_ Heated Purge: (Y/N) N  
 Calibration Start Date: 10/10/2017 09:35 Calibration End Date: 10/10/2017 09:58 Calibration ID: 30558

Calibration Files:

LEVEL:	LAB SAMPLE ID:	LAB FILE ID:
Level 1	STD001 280-390728/3	hfpo717J10026.d
Level 2	STD002 280-390728/4	hfpo717J10027.d
Level 3	STD003 280-390728/5	hfpo717J10028.d
Level 4	STD004 280-390728/6	hfpo717J10029.d
Level 5	STD005 280-390728/7	hfpo717J10030.d
Level 6	STD006 280-390728/8	hfpo717J10031.d
Level 7	STD007 280-390728/9	hfpo717J10032.d
Level 8	STD008 280-390728/10	hfpo717J10033.d

ANALYTE	CURVE TYPE	RESPONSE					CONCENTRATION (UG/L)				
		LVL 1 LVL 6	LVL 2 LVL 7	LVL 3 LVL 8	LVL 4	LVL 5	LVL 1 LVL 6	LVL 2 LVL 7	LVL 3 LVL 8	LVL 4	LVL 5
13C3 HFPO-DA	Ave	730749	745227	750427	718028	744600	10.0	10.0	10.0	10.0	10.0
		731935	729188	701420			10.0	10.0	10.0		

Curve Type Legend:

Ave = Average



FORM VI  
LCMS BY INTERNAL STANDARD - INITIAL CALIBRATION DATA  
RESPONSE AND CONCENTRATION

Lab Name: TestAmerica Denver Job No.: 280-105698-1 Analy Batch No.: 390728  
SDG No.: \_\_\_\_\_  
Instrument ID: LC\_LCMS7 GC Column: Synergi Hyd ID: \_\_\_\_\_ Heated Purge: (Y/N) N  
Calibration Start Date: 10/10/2017 09:35 Calibration End Date: 10/10/2017 09:58 Calibration ID: 30558

Calibration Files:

LEVEL:	LAB SAMPLE ID:	LAB FILE ID:
Level 1	STD001 280-390728/3	hfpo717J10026.d
Level 2	STD002 280-390728/4	hfpo717J10027.d
Level 3	STD003 280-390728/5	hfpo717J10028.d
Level 4	STD004 280-390728/6	hfpo717J10029.d
Level 5	STD005 280-390728/7	hfpo717J10030.d
Level 6	STD006 280-390728/8	hfpo717J10031.d
Level 7	STD007 280-390728/9	hfpo717J10032.d
Level 8	STD008 280-390728/10	hfpo717J10033.d

ANALYTE	IS REF	CURVE TYPE	RESPONSE					CONCENTRATION (UG/L)				
			LVL 1 LVL 6	LVL 2 LVL 7	LVL 3 LVL 8	LVL 4	LVL 5	LVL 1 LVL 6	LVL 2 LVL 7	LVL 3 LVL 8	LVL 4	LVL 5
HFPO-DA	13CP ODA	Lin1	31020 739399	63823 1790812	89272 3654104	167109	378047	0.250 10.0	0.500 25.0	1.00 50.0	2.00	5.00

Curve Type Legend:

[Lin1 = Linear 1/conc ISTD

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10026.d  
Lims ID: std001  
Client ID:  
Sample Type: IC Calib Level: 1  
Inject. Date: 10-Oct-2017 09:35:28 ALS Bottle#: 2 Worklist Smp#: 3  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: L1  
Misc. Info.: HFPO17J10  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Sublist: chrom-HFPO\*sub1  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 10-Oct-2017 12:51:45 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM  
Process Host: XAWRK005

First Level Reviewer: meyer Date: 10-Oct-2017 11:50:42

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 > 286.8	0.880	0.880	0.0	1.000	730749	10.0	397	
---------------	-------	-------	-----	-------	--------	------	-----	--

\* 2 13C3 HFPO-DA (IS)

331.8 > 286.8	0.880	0.880	0.0		730749	10.0	397	
---------------	-------	-------	-----	--	--------	------	-----	--

1 Perfluoro(2-propoxypropanoic) acid

328.8 > 284.8	0.893	0.885	0.008	1.000	31020	0.2036	14.1	M
---------------	-------	-------	-------	-------	-------	--------	------	---

## QC Flag Legend

Review Flags

M - Manually Integrated

## Reagents:

HFPO\_CAL-1\_00031

Amount Added: 1.00

Units: mL

TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10026.d

Injection Date: 10-Oct-2017 09:35:28

Instrument ID: LC\_LCMS7

Lims ID: std001

Client ID:

Operator ID: JBH

ALS Bottle#:

2

Worklist Smp#:

3

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid (M)

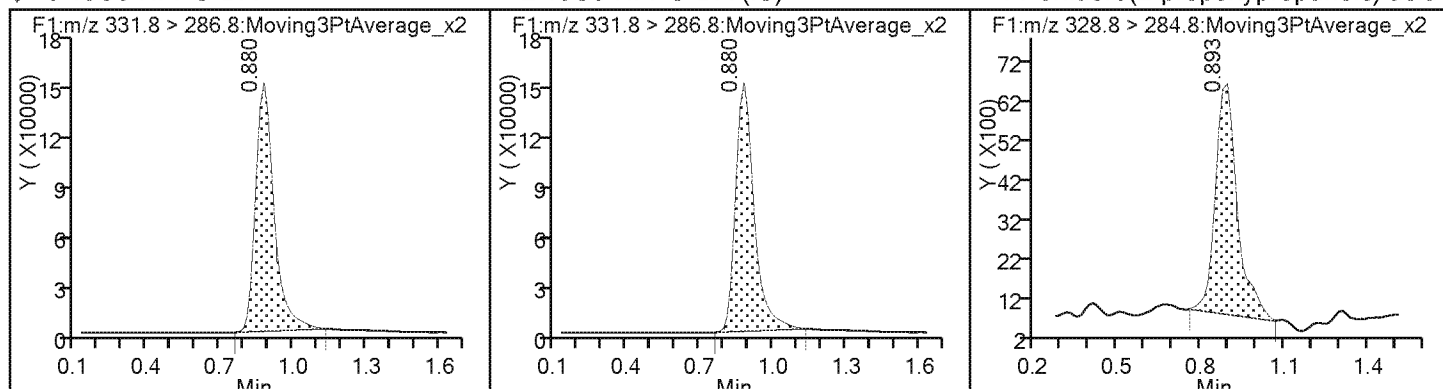


## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfp0717J10026.d  
Injection Date: 10-Oct-2017 09:35:28 Instrument ID: LC\_LCMS7  
Lims ID: std001  
Client ID:  
Operator ID: JBH ALS Bottle#: 2 Worklist Smp#: 3  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Method: HFPO Limit Group: LC - 8321A\_HFPO\_Du  
Column: Detector F1:MRM

1 Perfluoro(2-propoxypropanoic) acid, CAS: 13252-13-6  
Signal: 1


RT: 0.89  
Area: 24407  
Amount: 0.162386  
Amount Units: ug/l

## Processing Integration Results

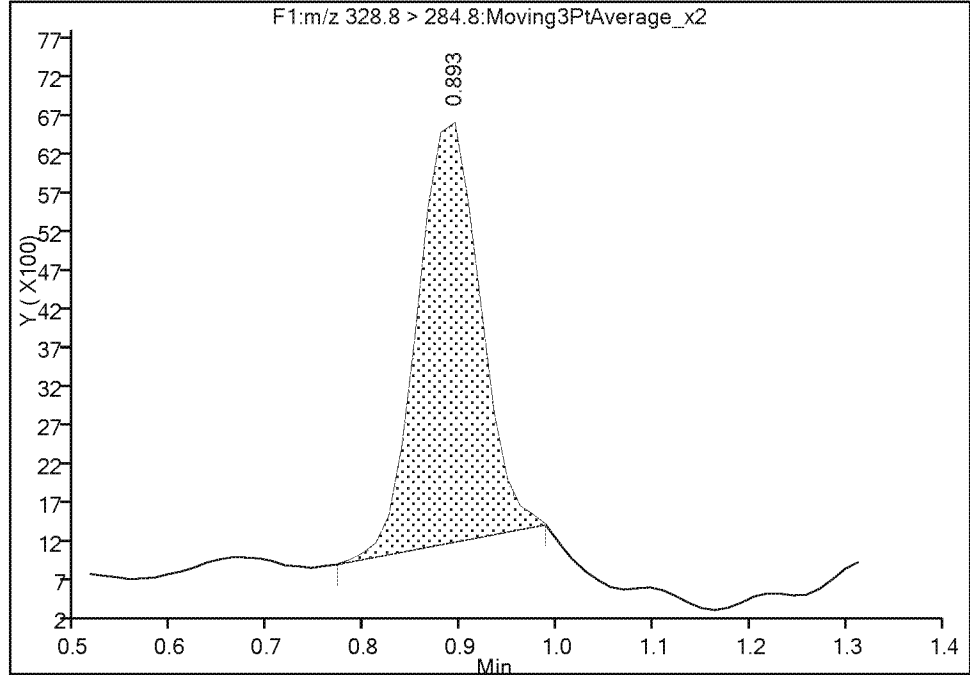

RT: 0.89  
Area: 31020  
Amount: 0.203553  
Amount Units: ug/l

## Manual Integration Results

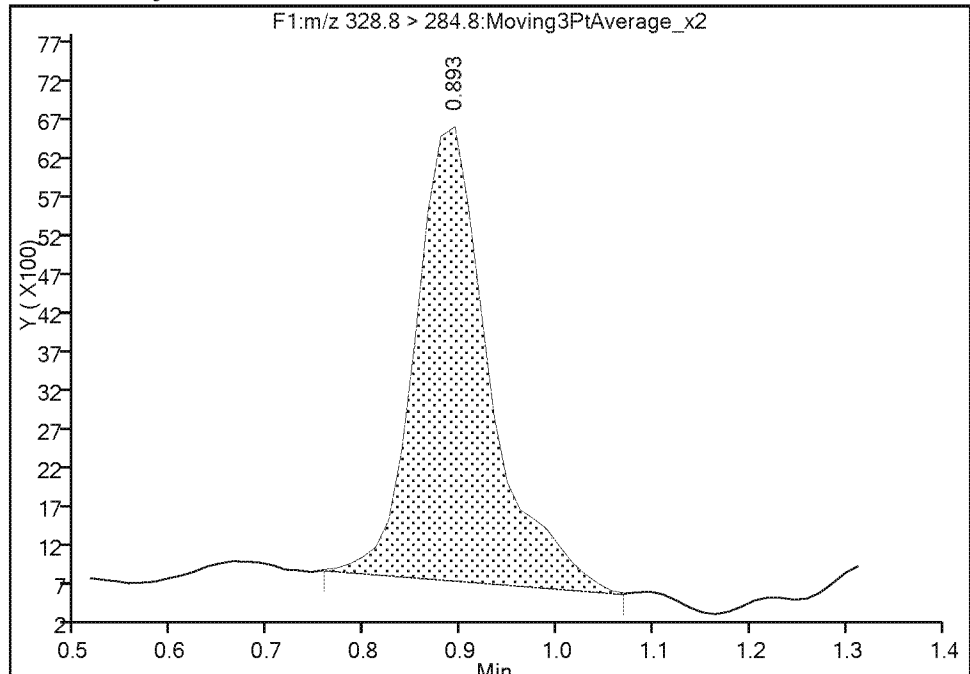

Reviewer: meyer, 10-Oct-2017 11:50:40  
Audit Action: Manually Integrated

Audit Reason: Baseline  
Page 152 of 237

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10027.d  
Lims ID: std002  
Client ID:  
Sample Type: IC Calib Level: 2  
Inject. Date: 10-Oct-2017 09:38:42 ALS Bottle#: 3 Worklist Smp#: 4  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: L2  
Misc. Info.: HFPO17J10  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Sublist: chrom-HFPO\*sub1  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 10-Oct-2017 12:51:46 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM  
Process Host: XAWRK005

First Level Reviewer: meyer Date: 10-Oct-2017 11:50:49

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\* 2 13C3 HFPO-DA (IS)

331.8 > 286.8 0.880 0.880 0.0 745227 10.0 452

\$ 3 13C3 HFPO-DA

331.8 > 286.8 0.880 0.880 0.0 1.000 745227 10.2 452

1 Perfluoro(2-propoxypropanoic) acid

328.8 > 284.8 0.880 0.885 -0.005 1.000 63823 0.6303 36.5

**Reagents:**

HFPO\_CAL-2\_00032

Amount Added: 1.00

Units: mL

TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10027.d

Injection Date: 10-Oct-2017 09:38:42

Instrument ID: LC\_LCMS7

Lims ID: std002

Client ID:

Operator ID: JBH

ALS Bottle#: 3

Worklist Smp#: 4

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\* 2 13C3 HFPO-DA (IS)

\$ 3 13C3 HFPO-DA

1 Perfluoro(2-propoxypropanoic) acid

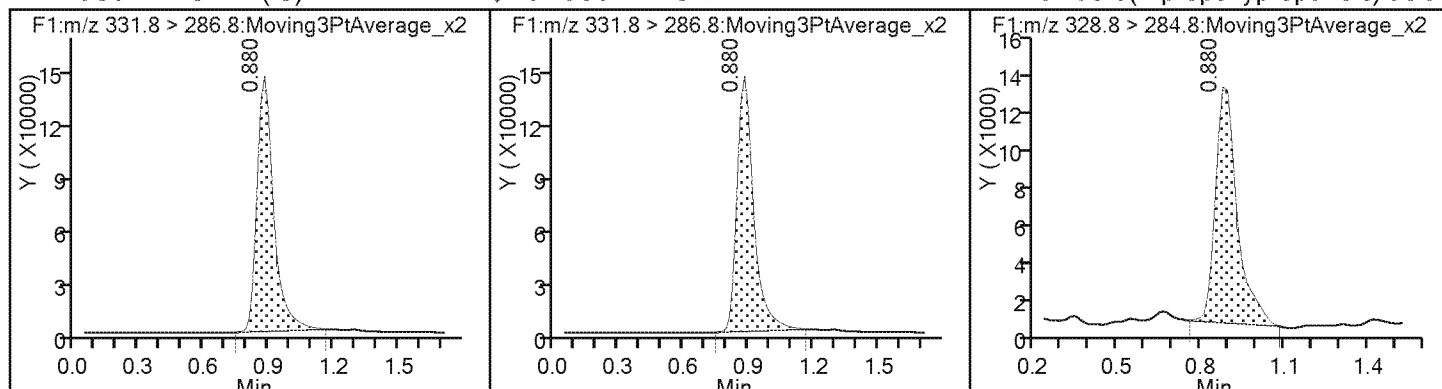

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10028.d  
Lims ID: std003  
Client ID:  
Sample Type: IC Calib Level: 3  
Inject. Date: 10-Oct-2017 09:41:56 ALS Bottle#: 4 Worklist Smp#: 5  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: L3  
Misc. Info.: HFPO17J10  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Sublist: chrom-HFPO\*sub1  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 10-Oct-2017 12:51:47 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM  
Process Host: XAWRK005

First Level Reviewer: meyer Date: 10-Oct-2017 11:50:52

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 > 286.8	0.880	0.880	0.0	1.000	750427	10.3	417	
---------------	-------	-------	-----	-------	--------	------	-----	--

\* 2 13C3 HFPO-DA (IS)

331.8 > 286.8	0.880	0.880	0.0		750427	10.0	417	
---------------	-------	-------	-----	--	--------	------	-----	--

1 Perfluoro(2-propoxypropanoic) acid

328.8 > 284.8	0.880	0.885	-0.005	1.000	89272	0.9595	50.3	
---------------	-------	-------	--------	-------	-------	--------	------	--

**Reagents:**

HFPO\_CAL-3\_00031 Amount Added: 1.00 Units: mL

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10028.d

Injection Date: 10-Oct-2017 09:41:56

Instrument ID: LC\_LCMS7

Lims ID: std003

Client ID:

Operator ID: JBH

ALS Bottle#:

4

Worklist Smp#:

5

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

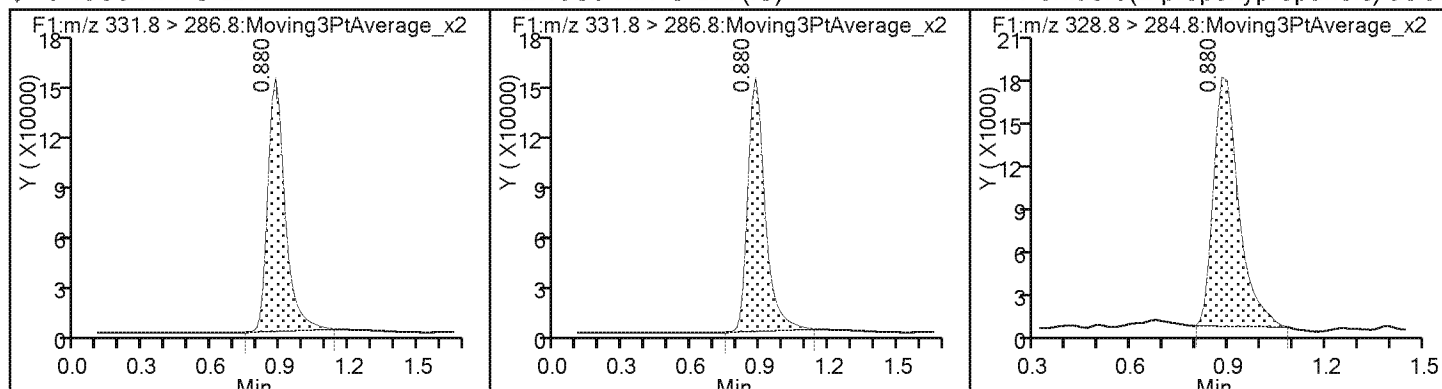



TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10029.d  
Lims ID: std004  
Client ID:  
Sample Type: IC Calib Level: 4  
Inject. Date: 10-Oct-2017 09:45:11 ALS Bottle#: 5 Worklist Smp#: 6  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: L4  
Misc. Info.: HFPO17J10  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Sublist: chrom-HFPO\*sub1  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 10-Oct-2017 12:51:47 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM  
Process Host: XAWRK005

First Level Reviewer: meyera Date: 10-Oct-2017 11:50:55

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\* 2 13C3 HFPO-DA (IS)

331.8 > 286.8 0.880 0.880 0.0 718028 10.0 438

\$ 3 13C3 HFPO-DA

331.8 > 286.8 0.880 0.880 0.0 1.000 718028 9.82 438

1 Perfluoro(2-propoxypropanoic) acid

328.8 > 284.8 0.880 0.885 -0.005 1.000 167109 2.08 143

**Reagents:**

HFPO\_CAL-4\_00031 Amount Added: 1.00 Units: mL

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10029.d

Injection Date: 10-Oct-2017 09:45:11

Instrument ID: LC\_LCMS7

Lims ID: std004

Client ID:

Operator ID: JBH

ALS Bottle#:

5

Worklist Smp#:

6

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\* 2 13C3 HFPO-DA (IS)

\$ 3 13C3 HFPO-DA

1 Perfluoro(2-propoxypropanoic) acid

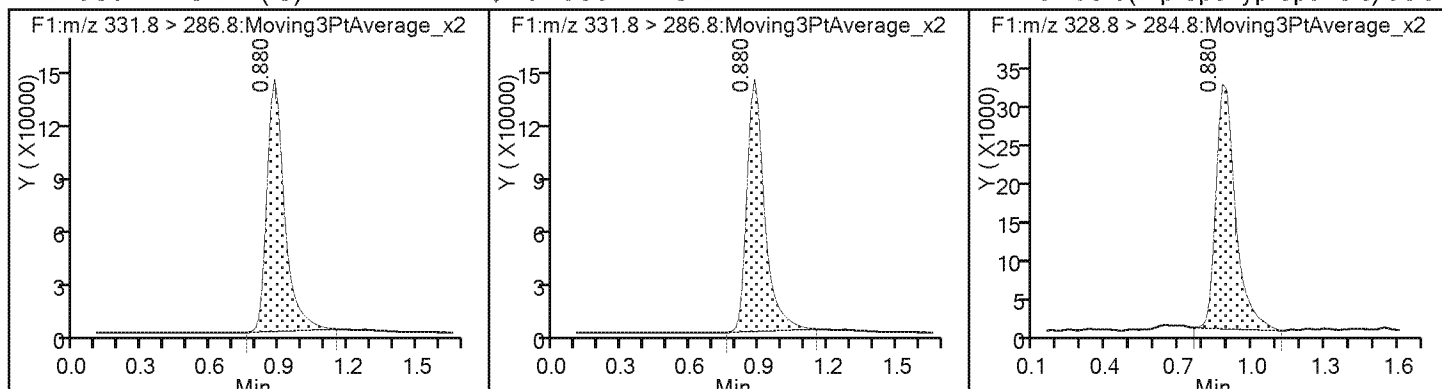

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10030.d  
Lims ID: std005  
Client ID:  
Sample Type: IC Calib Level: 5  
Inject. Date: 10-Oct-2017 09:48:25 ALS Bottle#: 6 Worklist Smp#: 7  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: L5  
Misc. Info.: HFPO17J10  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Sublist: chrom-HFPO\*sub1  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 10-Oct-2017 12:51:48 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM  
Process Host: XAWRK005

First Level Reviewer: meyer Date: 10-Oct-2017 11:50:57

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 > 286.8	0.880	0.880	0.0	1.000	744600	10.2	433	
---------------	-------	-------	-----	-------	--------	------	-----	--

\* 2 13C3 HFPO-DA (IS)

331.8 > 286.8	0.880	0.880	0.0		744600	10.0	433	
---------------	-------	-------	-----	--	--------	------	-----	--

1 Perfluoro(2-propoxypropanoic) acid

328.8 > 284.8	0.893	0.885	0.008	1.000	378047	4.80	223	
---------------	-------	-------	-------	-------	--------	------	-----	--

**Reagents:**

HFPO\_CAL-5\_00070

Amount Added: 1.00

Units: mL

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10030.d

Injection Date: 10-Oct-2017 09:48:25

Instrument ID: LC\_LCMS7

Lims ID: std005

Client ID:

Operator ID: JBH

ALS Bottle#:

6

Worklist Smp#:

7

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

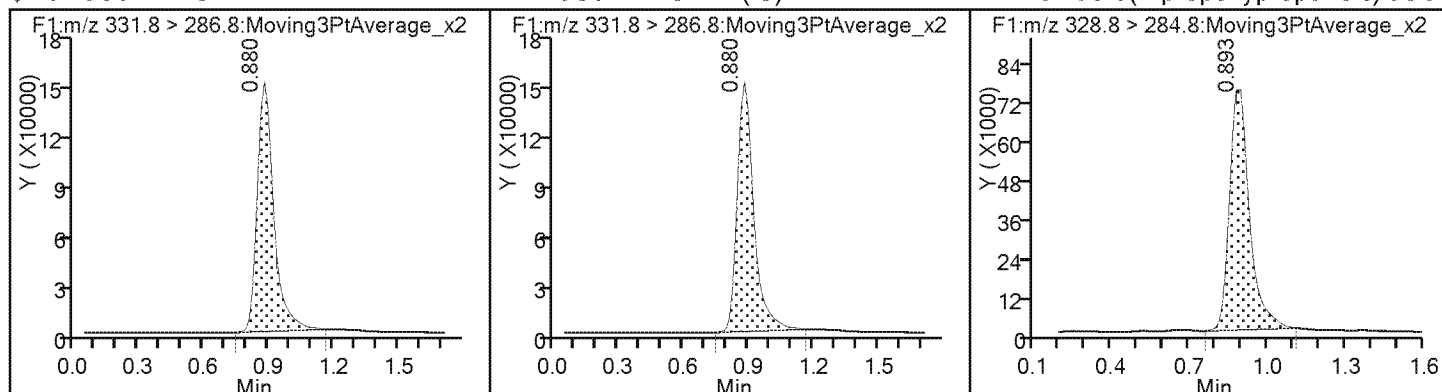

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10031.d  
Lims ID: std006  
Client ID:  
Sample Type: IC Calib Level: 6  
Inject. Date: 10-Oct-2017 09:51:39 ALS Bottle#: 7 Worklist Smp#: 8  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: L6  
Misc. Info.: HFPO17J10  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Sublist: chrom-HFPO\*sub1  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 10-Oct-2017 12:51:49 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM  
Process Host: XAWRK005

First Level Reviewer: meyer Date: 10-Oct-2017 11:51:00

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\* 2 13C3 HFPO-DA (IS)

331.8 > 286.8 0.880 0.880 0.0 731935 10.0 379

\$ 3 13C3 HFPO-DA

331.8 > 286.8 0.880 0.880 0.0 1.000 731935 10.0 379

1 Perfluoro(2-propoxypropanoic) acid

328.8 > 284.8 0.880 0.885 -0.005 1.000 739399 9.77 298

**Reagents:**

HFPO\_CAL-6\_00070

Amount Added: 1.00

Units: mL

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10031.d

Injection Date: 10-Oct-2017 09:51:39

Instrument ID: LC\_LCMS7

Lims ID: std006

Client ID:

Operator ID: JBH

ALS Bottle#: 7

Worklist Smp#: 8

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\* 2 13C3 HFPO-DA (IS)

\$ 3 13C3 HFPO-DA

1 Perfluoro(2-propoxypropanoic) acid

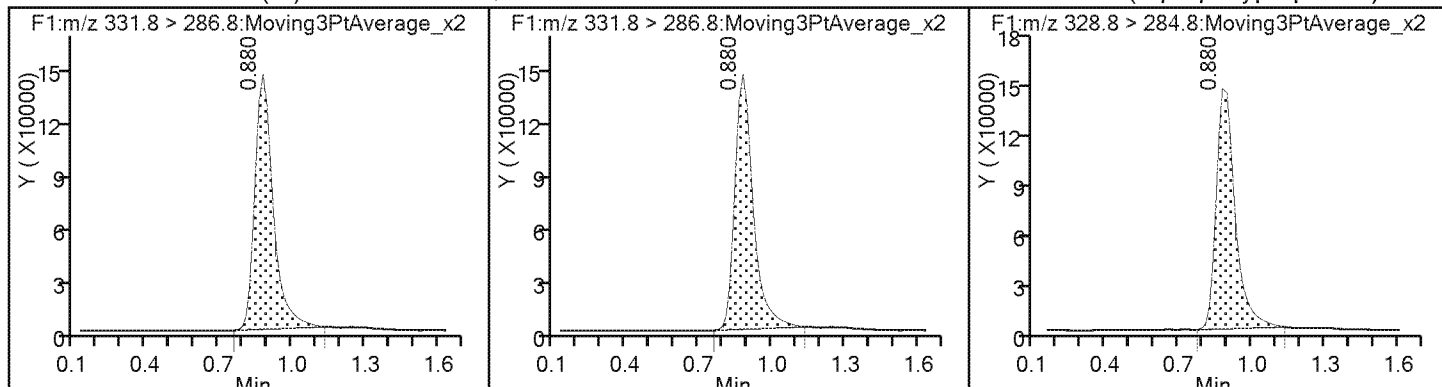

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10032.d  
Lims ID: std007  
Client ID:  
Sample Type: IC Calib Level: 7  
Inject. Date: 10-Oct-2017 09:54:53 ALS Bottle#: 8 Worklist Smp#: 9  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: L7  
Misc. Info.: HFPO17J10  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Sublist: chrom-HFPO\*sub1  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 10-Oct-2017 12:51:50 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM  
Process Host: XAWRK005

First Level Reviewer: meyer Date: 10-Oct-2017 11:51:04

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 > 286.8	0.880	0.880	0.0	1.000	729188	9.97	404	
---------------	-------	-------	-----	-------	--------	------	-----	--

\* 2 13C3 HFPO-DA (IS)

331.8 > 286.8	0.880	0.880	0.0		729188	10.0	404	
---------------	-------	-------	-----	--	--------	------	-----	--

1 Perfluoro(2-propoxypropanoic) acid

328.8 > 284.8	0.880	0.885	-0.005	1.000	1790812	24.0	386	
---------------	-------	-------	--------	-------	---------	------	-----	--

**Reagents:**

HFPO\_CAL-7\_00031 Amount Added: 1.00 Units: mL

TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10032.d

Injection Date: 10-Oct-2017 09:54:53

Instrument ID: LC\_LCMS7

Lims ID: std007

Client ID:

Operator ID: JBH

ALS Bottle#:

8

Worklist Smp#:

9

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

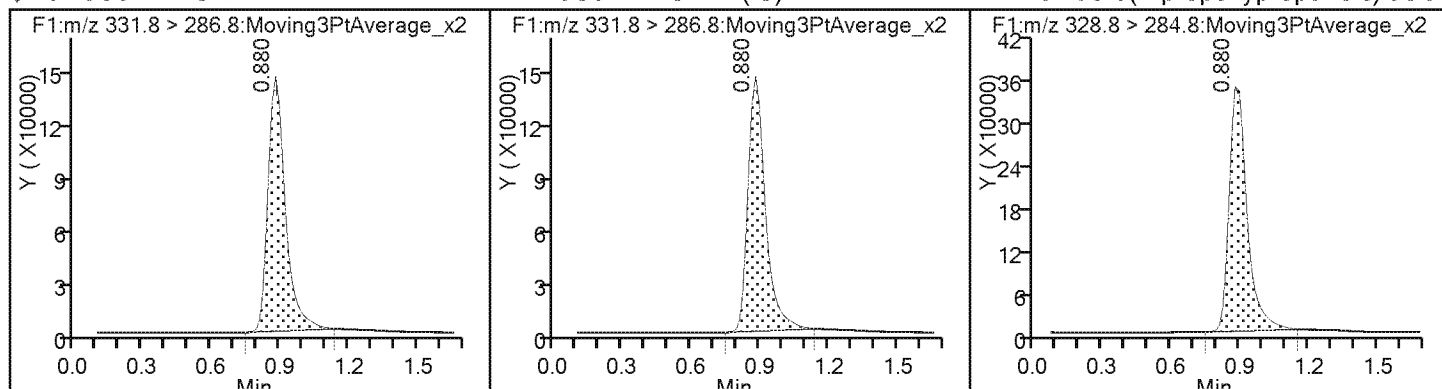



TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Lims ID: std008  
Client ID:  
Sample Type: IC Calib Level: 8  
Inject. Date: 10-Oct-2017 09:58:07 ALS Bottle#: 9 Worklist Smp#: 10  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: L8  
Misc. Info.: HFPO17J10  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Sublist: chrom-HFPO\*sub1  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 10-Oct-2017 12:51:51 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM  
Process Host: XAWRK005

First Level Reviewer: meyer Date: 10-Oct-2017 11:51:08

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\* 2 13C3 HFPO-DA (IS)

331.8 > 286.8 0.880 0.880 0.0 701420 10.0 373

\$ 3 13C3 HFPO-DA

331.8 > 286.8 0.880 0.880 0.0 1.000 701420 9.59 373

1 Perfluoro(2-propoxypropanoic) acid

328.8 > 284.8 0.893 0.885 0.008 1.000 3654104 51.3 421

**Reagents:**

HFPO\_CAL-8\_00031 Amount Added: 1.00 Units: mL

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d

Injection Date: 10-Oct-2017 09:58:07

Instrument ID: LC\_LCMS7

Lims ID: std008

Client ID:

Operator ID: JBH

ALS Bottle#:

9

Worklist Smp#:

10

Injection Vol: 20.0 ul

Dil. Factor:

1.0000

Method: HFPO


Limit Group:

LC - 8321A\_HFPO\_Du

\* 2 13C3 HFPO-DA (IS)

\$ 3 13C3 HFPO-DA

1 Perfluoro(2-propoxypropanoic) acid

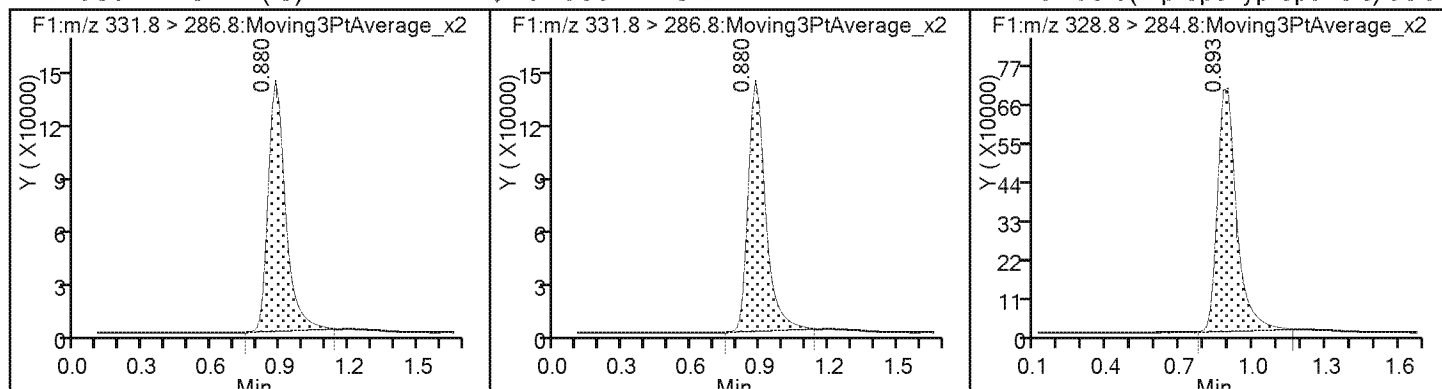

FORM VII  
LCMS CONTINUING CALIBRATION DATA

Lab Name: TestAmerica Denver Job No.: 280-105698-1  
SDG No.: \_\_\_\_\_  
Lab Sample ID: ICV 280-387775/13 Calibration Date: 09/14/2017 15:10  
Instrument ID: LC\_LCMS7 Calib Start Date: 09/14/2017 14:40  
GC Column: Synergi Hydro ID: \_\_\_\_\_ Calib End Date: 09/14/2017 15:01  
Lab File ID: hfpo717I14062.d Conc. Units: ug/L

ANALYTE	CURVE TYPE	AVE RRF	RRF	MIN RRF	CALC AMOUNT	SPIKE AMOUNT	%D	MAX %D
Perfluoro(2-propoxypropanoic ) acid	Lin1		0.9462		1.89	2.00	-5.3	20.0
13C3 HFPO-DA	Ave	192740	197806		10.3	10.0	2.6	

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\hfpo717114062.d  
Lims ID: ICV  
Client ID:  
Sample Type: ICV  
Inject. Date: 14-Sep-2017 15:10:31 ALS Bottle#: 10 Worklist Smp#: 13  
Injection Vol: 10.0 ul Dil. Factor: 1.0000  
Sample Info: ICV  
Misc. Info.: HFPO17114  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Sublist:

Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 15-Sep-2017 07:29:44 Calib Date: 14-Sep-2017 15:01:22  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\hfpo717114059.d

Column 1 : Det: F1:MRM  
Process Host: XAWRK034

First Level Reviewer: meyera Date: 15-Sep-2017 07:28:43

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 > 286.8 0.988 0.981 0.007 1.000 1978058 10.3 436

\* 2 13C3 HFPO-DA (IS)

331.8 > 286.8 0.988 0.981 0.007 1978058 10.0 436

1 Perfluoro(2-propoxypropanoic) acid

328.8 > 284.8 0.988 0.986 0.002 1.000 374307 1.89 162

**Reagents:**

HFPO\_ICV\_00031 Amount Added: 1.00 Units: mL

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20170915-62647.b\hfpo717114062.d

Injection Date: 14-Sep-2017 15:10:31

Instrument ID: LC\_LCMS7

Lims ID: ICV

Client ID:

Operator ID: JBH

ALS Bottle#: 10

Worklist Smp#: 13

Injection Vol: 10.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

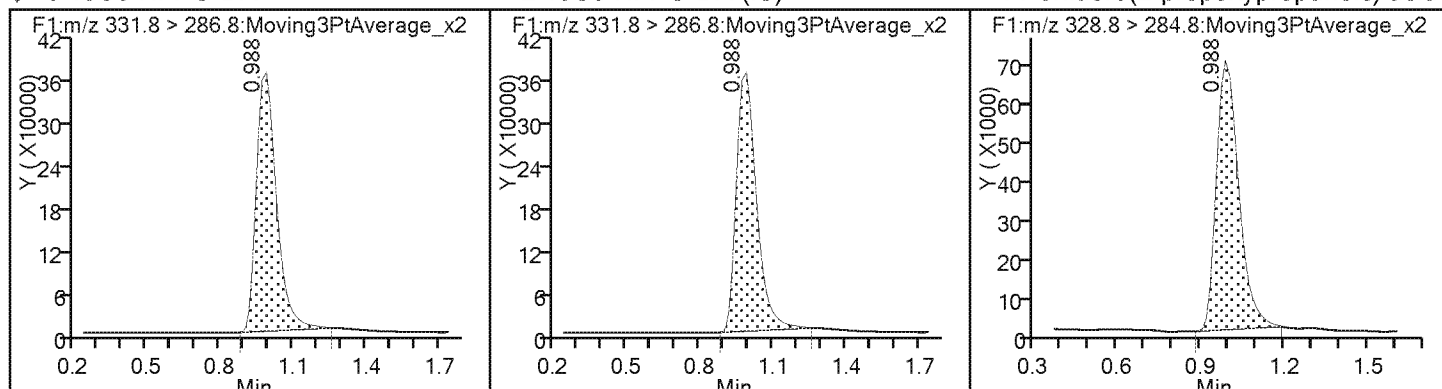

FORM VII  
LCMS CONTINUING CALIBRATION DATA

Lab Name: TestAmerica Denver Job No.: 280-105698-1  
SDG No.: \_\_\_\_\_  
Lab Sample ID: ICV 280-390728/13 Calibration Date: 10/10/2017 10:07  
Instrument ID: LC\_LCMS7 Calib Start Date: 10/10/2017 09:35  
GC Column: Synergi Hydro ID: \_\_\_\_\_ Calib End Date: 10/10/2017 09:58  
Lab File ID: hfpo717J10036.d Conc. Units: ug/L

ANALYTE	CURVE TYPE	AVE RRF	RRF	MIN RRF	CALC AMOUNT	SPIKE AMOUNT	%D	MAX %D
HFPO-DA	Lin1		1.154		2.07	2.00	3.3	20.0
13C3 HFPO-DA	Ave	73145	72923		9.97	10.0	-0.3	

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10036.d  
Lims ID: ICV  
Client ID:  
Sample Type: ICV  
Inject. Date: 10-Oct-2017 10:07:48 ALS Bottle#: 10 Worklist Smp#: 13  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: ICV  
Misc. Info.: HFPO17J10  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Sublist:

Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 10-Oct-2017 12:51:53 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM  
Process Host: XAWRK005

First Level Reviewer: meyer Date: 10-Oct-2017 11:51:34

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 > 286.8	0.880	0.880	0.0	1.000	729225	9.97	396	
---------------	-------	-------	-----	-------	--------	------	-----	--

\* 2 13C3 HFPO-DA (IS)

331.8 > 286.8	0.880	0.880	0.0		729225	10.0	396	
---------------	-------	-------	-----	--	--------	------	-----	--

1 Perfluoro(2-propoxypropanoic) acid

328.8 > 284.8	0.893	0.885	0.008	1.000	168368	2.07	111	
---------------	-------	-------	-------	-------	--------	------	-----	--

**Reagents:**

HFPO\_ICV\_00032 Amount Added: 1.00 Units: mL

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10036.d

Injection Date: 10-Oct-2017 10:07:48

Instrument ID: LC\_LCMS7

Lims ID: ICV

Client ID:

Operator ID: JBH

ALS Bottle#: 10

Worklist Smp#: 13

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

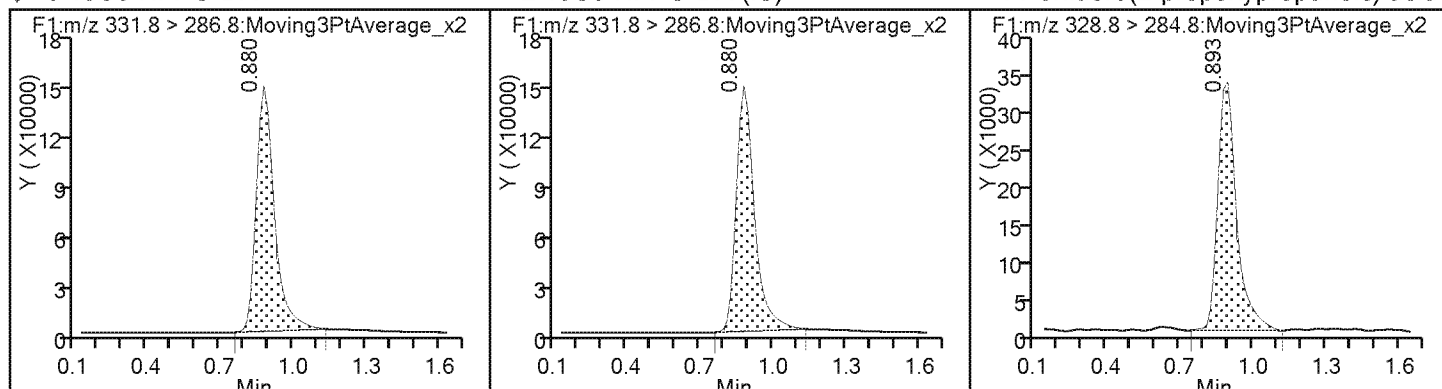



FORM VII  
LCMS CONTINUING CALIBRATION DATA

Lab Name: TestAmerica Denver Job No.: 280-105698-1  
SDG No.: \_\_\_\_\_  
Lab Sample ID: CCV 280-390728/24 Calibration Date: 10/10/2017 10:43  
Instrument ID: LC\_LCMS7 Calib Start Date: 10/10/2017 09:35  
GC Column: Synergi Hydro ID: \_\_\_\_\_ Calib End Date: 10/10/2017 09:58  
Lab File ID: hfpo717J10047.d Conc. Units: ug/L

ANALYTE	CURVE TYPE	AVE RRF	RRF	MIN RRF	CALC AMOUNT	SPIKE AMOUNT	%D	MAX %D
HFPO-DA	Lin1		1.012		9.78	10.0	-2.2	20.0
13C3 HFPO-DA	Ave	73145	68787		9.40	10.0	-6.0	

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10047.d  
Lims ID: CCV L6  
Client ID:  
Sample Type: CCV  
Inject. Date: 10-Oct-2017 10:43:29 ALS Bottle#: 7 Worklist Smp#: 24  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: CCV L6  
Misc. Info.: HFPO17J10  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Sublist: chrom-HFPO\*sub1  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 10-Oct-2017 12:52:02 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM  
Process Host: XAWRK005

First Level Reviewer: meyer Date: 10-Oct-2017 11:52:05

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 > 286.8 0.839 0.880 -0.041 1.000 687867 9.40 327

\* 2 13C3 HFPO-DA (IS)

331.8 > 286.8 0.839 0.880 -0.041 687867 10.0 327

1 Perfluoro(2-propoxypropanoic) acid

328.8 > 284.8 0.839 0.885 -0.046 1.000 696191 9.78 224

**Reagents:**

HFPO\_CAL-6\_00070

Amount Added: 1.00

Units: mL

TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10047.d

Injection Date: 10-Oct-2017 10:43:29

Instrument ID: LC\_LCMS7

Lims ID: CCV L6

Client ID:

Operator ID: JBH

ALS Bottle#: 7

Worklist Smp#: 24

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

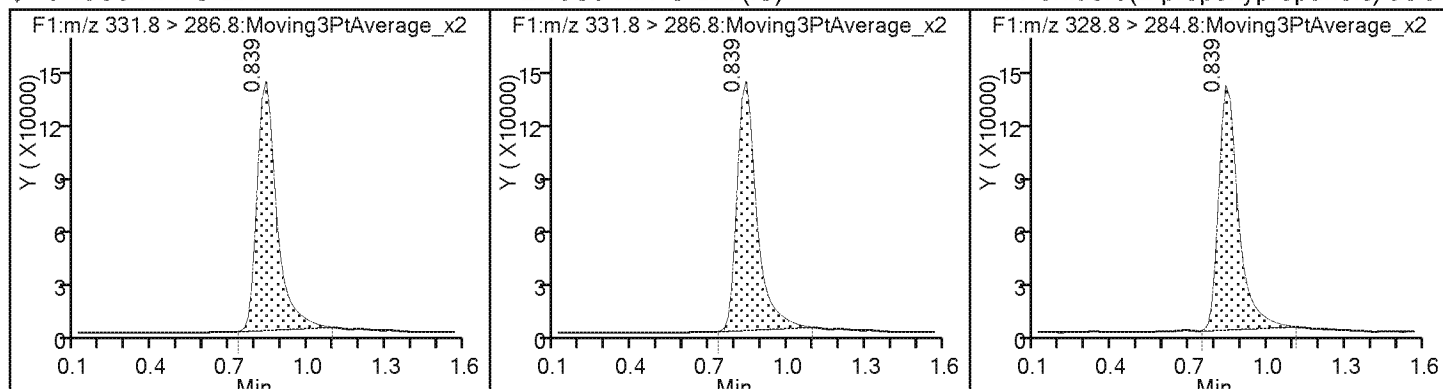

FORM VII  
LCMS CONTINUING CALIBRATION DATA

Lab Name: TestAmerica Denver Job No.: 280-105698-1  
SDG No.: \_\_\_\_\_  
Lab Sample ID: CCV 280-402806/3 Calibration Date: 01/25/2018 10:38  
Instrument ID: LC\_LCMS7 Calib Start Date: 10/10/2017 09:35  
GC Column: Synergi Hydro ID: \_\_\_\_\_ Calib End Date: 10/10/2017 09:58  
Lab File ID: hfpo718A25006.d Conc. Units: ug/L

ANALYTE	CURVE TYPE	AVE RRF	RRF	MIN RRF	CALC AMOUNT	SPIKE AMOUNT	%D	MAX %D
HFPO-DA	Lin1		0.9944		9.61	10.0	-3.9	20.0
13C3 HFPO-DA	Ave	73145	52305		7.15	10.0	-28.5	

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25006.d  
Lims ID: CCV L6  
Client ID:  
Sample Type: CCV  
Inject. Date: 25-Jan-2018 10:38:20 ALS Bottle#: 7 Worklist Smp#: 3  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: CCV L6  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Sublist: chrom-HFPO\*sub1  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:16 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyera Date: 25-Jan-2018 15:25:38

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 > 286.8 1.015 0.961 0.054 1.000 523048 7.15 1316

\* 2 13C3 HFPO-DA (IS)

331.8 > 286.8 1.015 0.961 0.054 523048 10.0 1316

1 Perfluoro(2-propoxypropanoic) acid

328.8 > 284.8 1.029 0.988 0.041 1.000 520124 9.61 87.5

**Reagents:**

HFPO\_CAL-6\_00078

Amount Added: 1.00

Units: mL

TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25006.d

Injection Date: 25-Jan-2018 10:38:20

Instrument ID: LC\_LCMS7

Lims ID: CCV L6

Client ID:

Operator ID: JBH

ALS Bottle#:

7

Worklist Smp#:

3

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

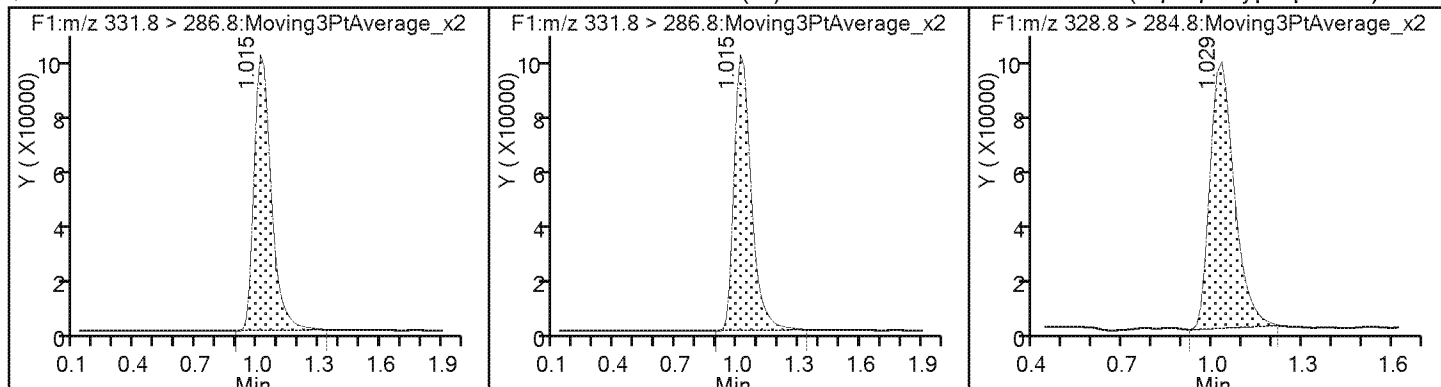

FORM VII  
LCMS CONTINUING CALIBRATION DATA

Lab Name: TestAmerica Denver Job No.: 280-105698-1  
SDG No.: \_\_\_\_\_  
Lab Sample ID: CCV 280-402806/11 Calibration Date: 01/25/2018 11:04  
Instrument ID: LC\_LCMS7 Calib Start Date: 10/10/2017 09:35  
GC Column: Synergi Hydro ID: \_\_\_\_\_ Calib End Date: 10/10/2017 09:58  
Lab File ID: hfpo718A25014.d Conc. Units: ug/L

ANALYTE	CURVE TYPE	AVE RRF	RRF	MIN RRF	CALC AMOUNT	SPIKE AMOUNT	%D	MAX %D
HFPO-DA	Lin1		1.058		5.01	5.00	0.2	20.0
13C3 HFPO-DA	Ave	73145	82916		11.3	10.0	13.4	

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25014.d  
Lims ID: CCV L5  
Client ID:  
Sample Type: CCV  
Inject. Date: 25-Jan-2018 11:04:18 ALS Bottle#: 6 Worklist Smp#: 11  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: CCV L5  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Sublist: chrom-HFPO\*sub1  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:23 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyera Date: 25-Jan-2018 15:26:39

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 > 286.8 0.961 0.961 0.0 1.000 829158 11.3 1627

\* 2 13C3 HFPO-DA (IS)

331.8 > 286.8 0.961 0.961 0.0 829158 10.0 1627

1 Perfluoro(2-propoxypropanoic) acid

328.8 > 284.8 0.961 0.988 -0.027 1.000 438620 5.01 149

**Reagents:**

HFPO\_CAL-5\_00078

Amount Added: 1.00

Units: mL



## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25014.d

Injection Date: 25-Jan-2018 11:04:18

Instrument ID: LC\_LCMS7

Lims ID: CCV L5

Client ID:

Operator ID: JBH

ALS Bottle#: 6

Worklist Smp#: 11

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

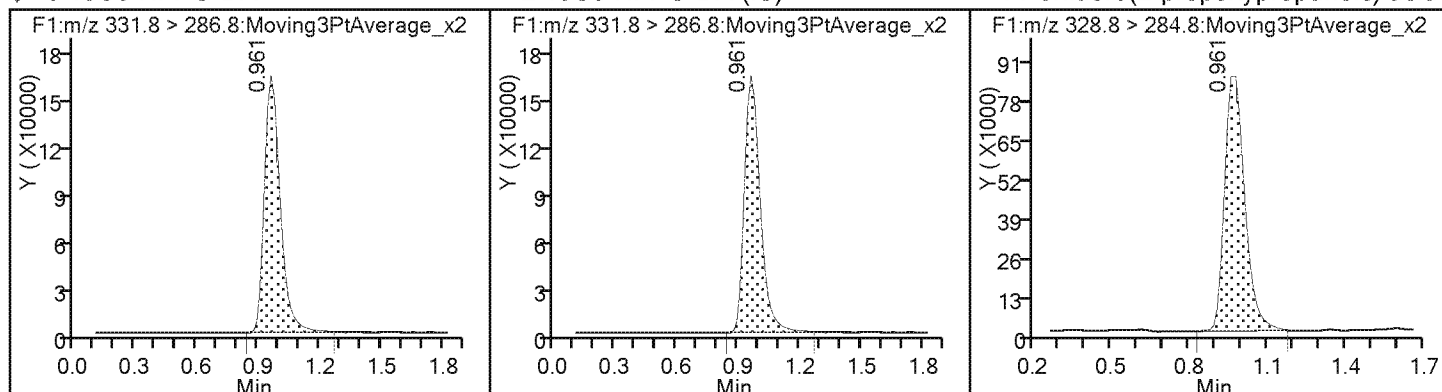

FORM VII  
LCMS CONTINUING CALIBRATION DATA

Lab Name: TestAmerica Denver Job No.: 280-105698-1  
SDG No.: \_\_\_\_\_  
Lab Sample ID: CCV 280-402806/19 Calibration Date: 01/25/2018 11:30  
Instrument ID: LC\_LCMS7 Calib Start Date: 10/10/2017 09:35  
GC Column: Synergi Hydro ID: \_\_\_\_\_ Calib End Date: 10/10/2017 09:58  
Lab File ID: hfpo718A25022.d Conc. Units: ug/L

ANALYTE	CURVE TYPE	AVE RRF	RRF	MIN RRF	CALC AMOUNT	SPIKE AMOUNT	%D	MAX %D
HFPO-DA	Lin1		1.020		9.86	10.0	-1.4	20.0
13C3 HFPO-DA	Ave	73145	78320		10.7	10.0	7.1	

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25022.d  
Lims ID: CCV L6  
Client ID:  
Sample Type: CCV  
Inject. Date: 25-Jan-2018 11:30:21 ALS Bottle#: 7 Worklist Smp#: 19  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: CCV L6  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Sublist: chrom-HFPO\*sub1  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:31 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyer Date: 25-Jan-2018 15:28:46

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 > 286.8 0.961 0.961 0.0 1.000 783195 10.7 1798

\* 2 13C3 HFPO-DA (IS)

331.8 > 286.8 0.961 0.961 0.0 783195 10.0 1798

1 Perfluoro(2-propoxypropanoic) acid

328.8 > 284.8 0.961 0.988 -0.027 1.000 799006 9.86 216

**Reagents:**

HFPO\_CAL-6\_00078

Amount Added: 1.00

Units: mL

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25022.d

Injection Date: 25-Jan-2018 11:30:21

Instrument ID: LC\_LCMS7

Lims ID: CCV L6

Client ID:

Operator ID: JBH

ALS Bottle#: 7

Worklist Smp#: 19

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

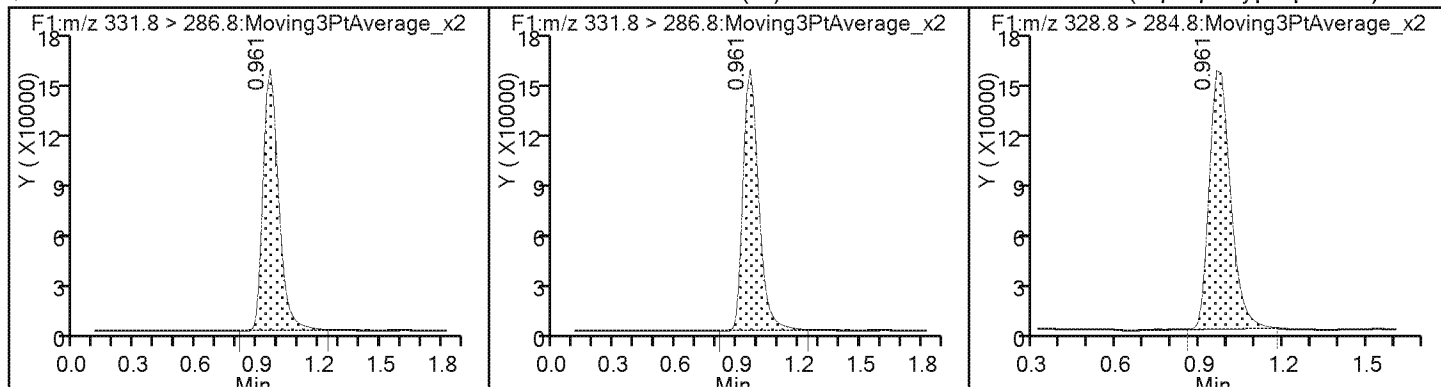

FORM VII  
LCMS CONTINUING CALIBRATION DATA

Lab Name: TestAmerica Denver Job No.: 280-105698-1  
SDG No.: \_\_\_\_\_  
Lab Sample ID: CCV 280-402806/28 Calibration Date: 01/25/2018 11:59  
Instrument ID: LC\_LCMS7 Calib Start Date: 10/10/2017 09:35  
GC Column: Synergi Hydro ID: \_\_\_\_\_ Calib End Date: 10/10/2017 09:58  
Lab File ID: hfpo718A25031.d Conc. Units: ug/L

ANALYTE	CURVE TYPE	AVE RRF	RRF	MIN RRF	CALC AMOUNT	SPIKE AMOUNT	%D	MAX %D
HFPO-DA	Lin1		1.079		5.11	5.00	2.3	20.0
13C3 HFPO-DA	Ave	73145	85353		11.7	10.0	16.7	

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25031.d  
Lims ID: CCV L5  
Client ID:  
Sample Type: CCV  
Inject. Date: 25-Jan-2018 11:59:50 ALS Bottle#: 6 Worklist Smp#: 28  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: CCV L5  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Sublist: chrom-HFPO\*sub1  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:39 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyera Date: 25-Jan-2018 15:29:07

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 > 286.8	0.961	0.961	0.0	1.000	853531	11.7	1158
---------------	-------	-------	-----	-------	--------	------	------

\* 2 13C3 HFPO-DA (IS)

331.8 > 286.8	0.961	0.961	0.0		853531	10.0	1158
---------------	-------	-------	-----	--	--------	------	------

1 Perfluoro(2-propoxypropanoic) acid

328.8 > 284.8	0.961	0.988	-0.027	1.000	460339	5.11	199
---------------	-------	-------	--------	-------	--------	------	-----

**Reagents:**

HFPO\_CAL-5\_00078

Amount Added: 1.00

Units: mL

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25031.d

Injection Date: 25-Jan-2018 11:59:50

Instrument ID: LC\_LCMS7

Lims ID: CCV L5

Client ID:

Operator ID: JBH

ALS Bottle#: 6

Worklist Smp#: 28

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)


1 Perfluoro(2-propoxypropanoic) acid

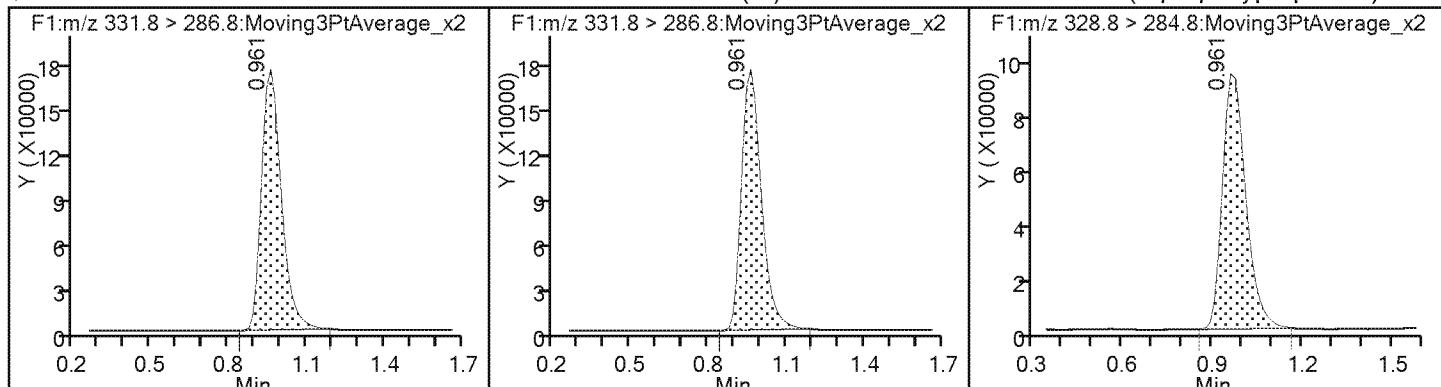

File: C:\MassLynx\6321.PRO\ACQUDB\HFPOMRM.lpr

Instrument: XEVO-TQMS\FBA453

Printed: Thursday, January 25, 2018 09:55:50 Mountain Standard Time

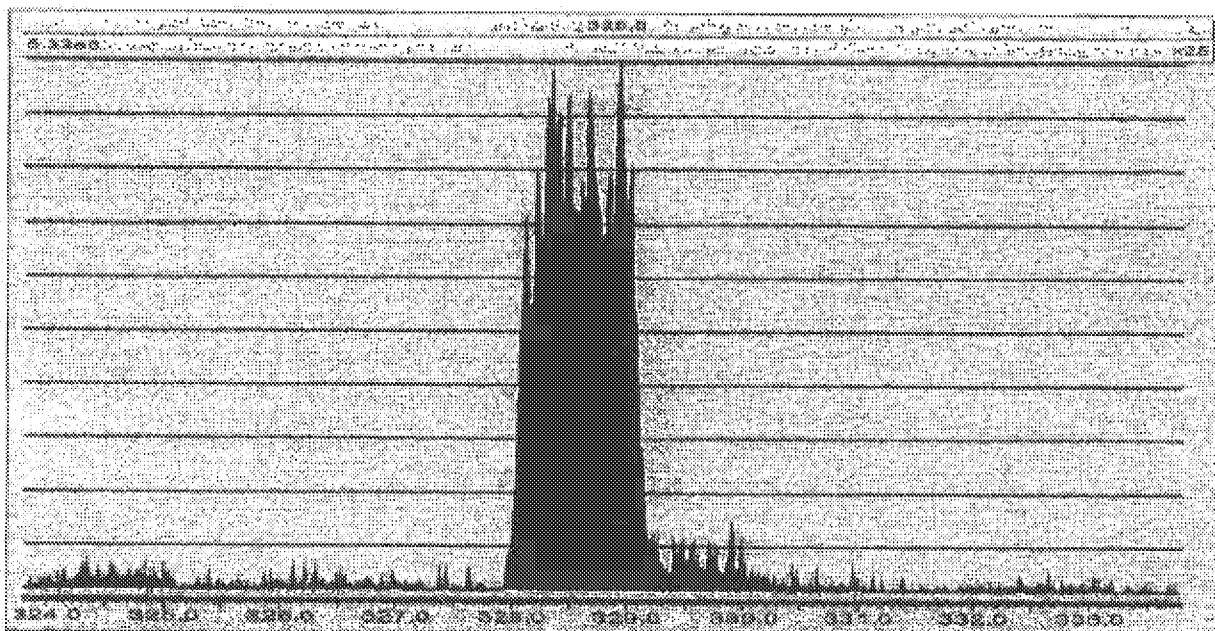

Type	Start Mass	End Mass	Set Mass
MS1 Scan	323.80	333.80	

Source (ES-)	Settings	Readbacks
Capillary (KV)	0.50	0.52
Cone (V)	10.00	-21.06
Extractor (V)	3.00	-10.61
Source Temperature (°C)	120	120
Desolvation Temperature (°C)	200	200
Cone Gas Flow (L/Hr)	50	50
Desolvation Gas Flow (L/Hr)	800	798
Collision Gas Flow (mL/Min)	0.15	0.03

Analyser	Settings	Readbacks
LM 1 Resolution	2.8	
HM 1 Resolution	14.8	
Ion Energy 1	0.7	
MS Mode Collision Energy	7.00	
MSMS Mode Collision Energy	20.00	
MS Mode Entrance	0.50	
MS Mode Exit	0.50	
Gas On MS Mode Entrance	0.50	
Gas On MS Mode Exit	0.50	
Gas On MSMS Mode Entrance	0.50	
Gas On MSMS Mode Exit	0.50	
Gas Off MS Mode Entrance	30.00	
Gas Off MS Mode Exit	30.00	
Gas Off MSMS Mode Entrance	2.00	
Gas Off MSMS Mode Exit	2.00	
ScanWave MS Mode Entrance	0.50	
ScanWave MS Mode Exit	0.50	
ScanWave MSMS Mode Entrance	0.50	
ScanWave MSMS Mode Exit	0.50	
LM 2 Resolution	2.9	
HM 2 Resolution	14.7	
Ion Energy 2	0.3	

Proton  
01/25/2018



File: C:\MassLynx\8321.PRO\ACQUDB\HFPOMRM.lpr

Instrument: XEVO-TQMS\FVBA453

Printed: Thursday, January 25, 2018 09:55:50 Mountain Standard Time

Multiplier 623.67  
Active Reservoir A

Pressure Gauges  
Collision Cell Pressure (mbar) 7.830201e-005

## Instrument Configuration

## Automatic Mode

MS Inter-scan delay (secs) 0.005  
Polarity/Mode switch Inter-scan delay (secs) 0.020  
Enhanced Inter-scan delay (secs) 0.020  
Inter-channel delay - See Tables

## MS 1 Delay Table:

R	delay
<= 0.500	0.005
<= 2.000	0.008
<= 4.000	0.010
<= 11.000	0.012
> 11.000	0.014

## MS 2 Delay Table:


R	delay
<= 8.000	0.005
<= 25.000	0.006
> 25.000	0.007

Phuong  
01/25/2018

File: C:\MassLynx\8321.PROVACQUDB\HFPOMRM.lpr

Instrument: XEVO-TQMS\FBA453

Printed: Thursday, January 25, 2018 09:56:24 Mountain Standard Time

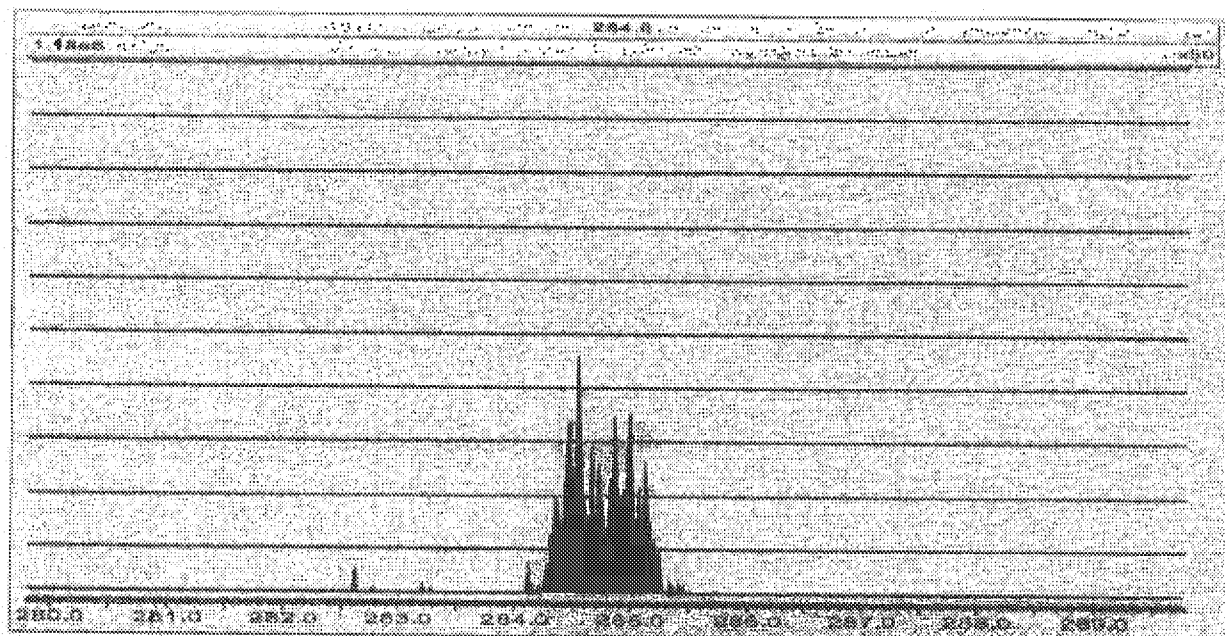

Type	Start Mass	End Mass	Set Mass
Daughter Scan	279.80	289.80	328.80

Source (ES-)	Settings	Readbacks
Capillary (kV)	0.50	0.52
Cone (V)	10.00	-21.08
Extractor (V)	3.00	-10.61
Source Temperature (°C)	120	120
Desolvation Temperature (°C)	200	200
Cone Gas Flow (L/Hr)	50	51
Desolvation Gas Flow (L/Hr)	800	793
Collision Gas Flow (mL/Min)	0.15	0.14

Analyser	Settings	Readbacks
LM 1 Resolution	2.8	
HM 1 Resolution	14.8	
Ion Energy 1	0.7	
MS Mode Collision Energy	7.00	
MSMS Mode Collision Energy	20.00	
MS Mode Entrance	0.50	
MS Mode Exit	0.50	
Gas On MS Mode Entrance	0.50	
Gas On MS Mode Exit	0.50	
Gas On MSMS Mode Entrance	0.50	
Gas On MSMS Mode Exit	0.50	
Gas Off MS Mode Entrance	30.00	
Gas Off MS Mode Exit	30.00	
Gas Off MSMS Mode Entrance	2.00	
Gas Off MSMS Mode Exit	2.00	
ScanWave MS Mode Entrance	0.50	
ScanWave MS Mode Exit	0.50	
ScanWave MSMS Mode Entrance	0.50	
ScanWave MSMS Mode Exit	0.50	
LM 2 Resolution	2.9	
HM 2 Resolution	14.7	
Ion Energy 2	0.3	

Phuriga  
4/14/2018

File: C:\MassLynx\8321.PROVACQUDB\HFPOMRM.lpr

Instrument: XEVO-TQMS\FBA453

Printed: Thursday, January 25, 2018 09:56:24 Mountain Standard Time

Multiplier 523.57  
Active Reservoir A

Pressure Gauges  
Collision Cell Pressure (mbar) 1.168540e-003

## Instrument Configuration

## Automatic Mode

MS Inter-scan delay (secs) 0.005

Polarity/Mode switch Inter-scan delay (secs) 0.020

Enhanced Inter-scan delay (secs) 0.020

Inter-channel delay - See Tables

## MS 1 Delay Table:

R	delay
<= 0.500	0.005
<= 2.000	0.006
<= 4.000	0.010
<= 11.000	0.012
> 11.000	0.014

## MS 2 Delay Table:

R	delay
<= 8.000	0.005
<= 25.000	0.006
> 25.000	0.007

Proton  
01/25/2018

## MS Method Report - MassLynx 4.1 SCN 843

Page 1 of 1

File: c:\masslynx\8321.pro\acq\db\hfpo.exp

Printed: Thursday, January 25, 2018 09:56:56 Mountain Standard Time

Creation Time Fri 18 Nov 2016 09:08:40  
Instrument Identifier XEVO-TQMS\FBA463  
Version Number 1.0  
Duration (min) 2.0  
Calibration Filename C:\MassLynx\IntelliStart\Results\Unit Mass Resolution\Calibration\_20100811

\_2.cal

Solvent Delay Divert Valve Enabled 0  
Number Of Functions 1

## Function 1 : MRM of 2 mass pairs, Time 0.00 to 2.00, ES-

Type MRM  
Ion Mode ES-  
Inter Channel Delay (sec) -1.000  
InterScan Time (sec) -1.000  
Span (Da) 0.5  
Start Time (min) 0.0  
End Time (min) 2.0

Ch	Prnt (Da)	Dau (Da)	Dwell (s)	Cons (V)	Coll (eV)	Delay (s)	Compound
1	328.80	284.80	0.400	10.00	7.00	-1.000	HFPO
2	331.80	286.80	0.400	10.00	7.00	-1.000	HFPO IS

Phurpa.  
01/25/2018

FORM I  
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: <u>TestAmerica Denver</u>	Job No.: <u>280-105698-1</u>
SDG No.: _____	
Client Sample ID: _____	Lab Sample ID: <u>MB 280-402648/1-A</u>
Matrix: <u>Water</u>	Lab File ID: <u>hfpo718A25007.d</u>
Analysis Method: <u>8321A</u>	Date Collected: _____
Extraction Method: <u>3535</u>	Date Extracted: <u>01/24/2018 15:20</u>
Sample wt/vol: <u>250 (mL)</u>	Date Analyzed: <u>01/25/2018 10:41</u>
Con. Extract Vol.: <u>5 (mL)</u>	Dilution Factor: <u>1</u>
Injection Volume: <u>20 (uL)</u>	GC Column: <u>Synergi Hydro</u> ID: _____
% Moisture: _____	GPC Cleanup: (Y/N) <u>N</u>
Analysis Batch No.: <u>402806</u>	Units: <u>ug/L</u>

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	<0.010		0.010	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	93		50-200

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25007.d  
Lims ID: MB 280-402648/1-A  
Client ID:  
Sample Type: MB  
Inject. Date: 25-Jan-2018 10:41:34 ALS Bottle#: 11 Worklist Smp#: 4  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: MB280-402648/1-A  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:16 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyera

Date: 25-Jan-2018 15:25:40

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 &gt; 286.8 1.015 0.961 0.054 1.000 677806 9.27 1613

\* 2 13C3 HFPO-DA (IS)

331.8 &gt; 286.8 1.015 0.961 0.054 677806 10.0 1613

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25007.d

Injection Date: 25-Jan-2018 10:41:34

Instrument ID: LC\_LCMS7

Lims ID: MB 280-402648/1-A

Client ID:

Operator ID: JBH

ALS Bottle#: 11

Worklist Smp#: 4

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid (ND)

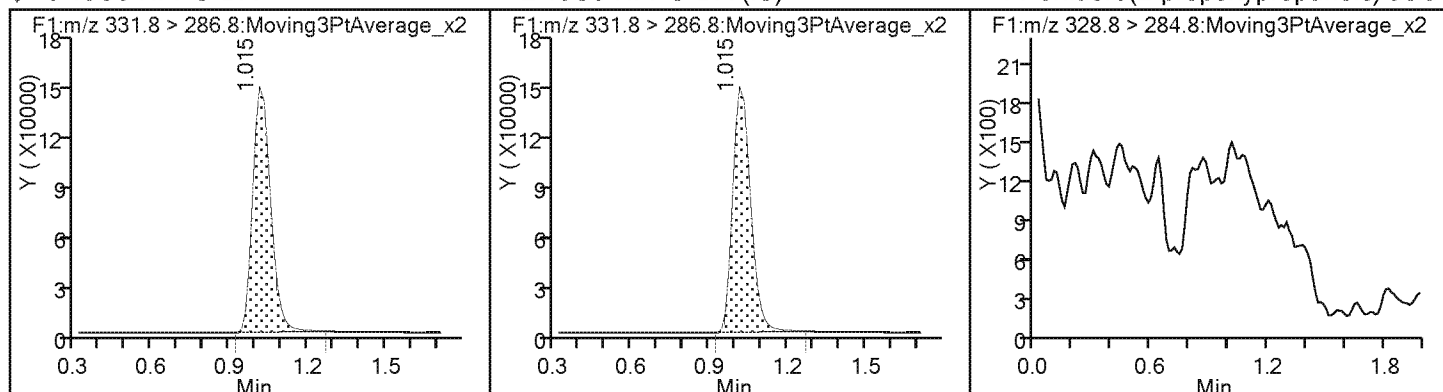

TestAmerica Denver  
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25007.d  
Lims ID: MB 280-402648/1-A  
Client ID:  
Sample Type: MB  
Inject. Date: 25-Jan-2018 10:41:34 ALS Bottle#: 11 Worklist Smp#: 4  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: MB280-402648/1-A  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:16 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyera

Date: 25-Jan-2018 15:25:40

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	9.27	92.67



FORM I  
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: TestAmerica Denver Job No.: 280-105698-1  
 SDG No.: \_\_\_\_\_  
 Client Sample ID: \_\_\_\_\_ Lab Sample ID: ICB 280-390728/11  
 Matrix: Water Lab File ID: hfpo717J10034.d  
 Analysis Method: 8321A Date Collected: \_\_\_\_\_  
 Extraction Method: \_\_\_\_\_ Date Extracted: \_\_\_\_\_  
 Sample wt/vol: 1(mL) Date Analyzed: 10/10/2017 10:01  
 Con. Extract Vol.: \_\_\_\_\_ Dilution Factor: 1  
 Injection Volume: 20(uL) GC Column: Synergi Hydro ID: \_\_\_\_\_  
 % Moisture: \_\_\_\_\_ GPC Cleanup: (Y/N) N  
 Analysis Batch No.: 390728 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	<0.50		0.50	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	100		50-200

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10034.d  
Lims ID: ICB  
Client ID:  
Sample Type: ICB  
Inject. Date: 10-Oct-2017 10:01:21 ALS Bottle#: 1 Worklist Smp#: 11  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: ICB  
Misc. Info.: HFPO17J10  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 10-Oct-2017 12:51:51 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK005

First Level Reviewer: meyer

Date: 10-Oct-2017 11:51:12

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 &gt; 286.8 0.880 0.880 0.0 1.000 732194 10.0 425

\* 2 13C3 HFPO-DA (IS)

331.8 &gt; 286.8 0.880 0.880 0.0 732194 10.0 425

1 Perfluoro(2-propoxypropanoic) acid

328.8 &gt; 284.8 0.880 0.885 -0.005 1.000 13993 -0.0270 8.1

## Reagents:

HFPO\_CAL-0\_00031

Amount Added: 1.00

Units: mL

TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10034.d

Injection Date: 10-Oct-2017 10:01:21

Instrument ID: LC\_LCMS7

Lims ID: ICB

Client ID:

Operator ID: JBH

ALS Bottle#: 1

Worklist Smp#: 11

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

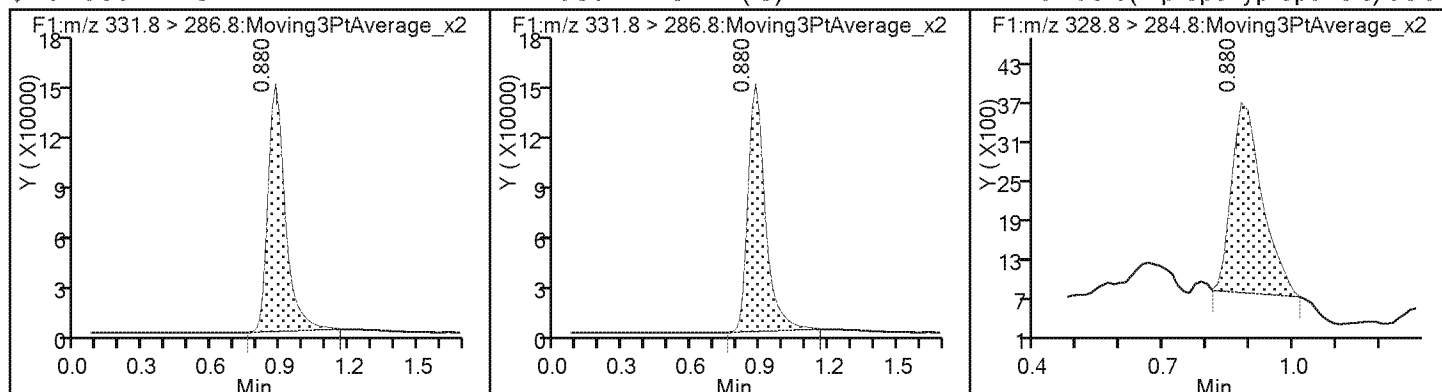

TestAmerica Denver  
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10034.d  
Lims ID: ICB  
Client ID:  
Sample Type: ICB  
Inject. Date: 10-Oct-2017 10:01:21 ALS Bottle#: 1 Worklist Smp#: 11  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: ICB  
Misc. Info.: HFPO17J10  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 10-Oct-2017 12:51:51 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK005

First Level Reviewer: meyer

Date: 10-Oct-2017 11:51:12

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	10.0	100.10

FORM I  
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: TestAmerica Denver Job No.: 280-105698-1  
 SDG No.: \_\_\_\_\_  
 Client Sample ID: \_\_\_\_\_ Lab Sample ID: LCS 280-402648/2-A  
 Matrix: Water Lab File ID: hfpo718A25008.d  
 Analysis Method: 8321A Date Collected: \_\_\_\_\_  
 Extraction Method: 3535 Date Extracted: 01/24/2018 15:20  
 Sample wt/vol: 250 (mL) Date Analyzed: 01/25/2018 10:44  
 Con. Extract Vol.: 5 (mL) Dilution Factor: 1  
 Injection Volume: 20 (uL) GC Column: Synergi Hydro ID: \_\_\_\_\_  
 % Moisture: \_\_\_\_\_ GPC Cleanup: (Y/N) N  
 Analysis Batch No.: 402806 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	0.193		0.010	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	91		50-200

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25008.d  
Lims ID: LCS 280-402648/2-A  
Client ID:  
Sample Type: LCS  
Inject. Date: 25-Jan-2018 10:44:48 ALS Bottle#: 12 Worklist Smp#: 5  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: LCS280-402648/2-A  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:16 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyera

Date: 25-Jan-2018 15:25:42

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 &gt; 286.8 0.961 0.961 0.0 1.000 665871 9.10 1913

\* 2 13C3 HFPO-DA (IS)

331.8 &gt; 286.8 0.961 0.961 0.0 665871 10.0 1913

1 Perfluoro(2-propoxypropanoic) acid

328.8 &gt; 284.8 0.975 0.988 -0.013 1.000 664038 9.64 420

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25008.d

Injection Date: 25-Jan-2018 10:44:48

Instrument ID: LC\_LCMS7

Lims ID: LCS 280-402648/2-A

Client ID:

Operator ID: JBH

ALS Bottle#: 12

Worklist Smp#: 5

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

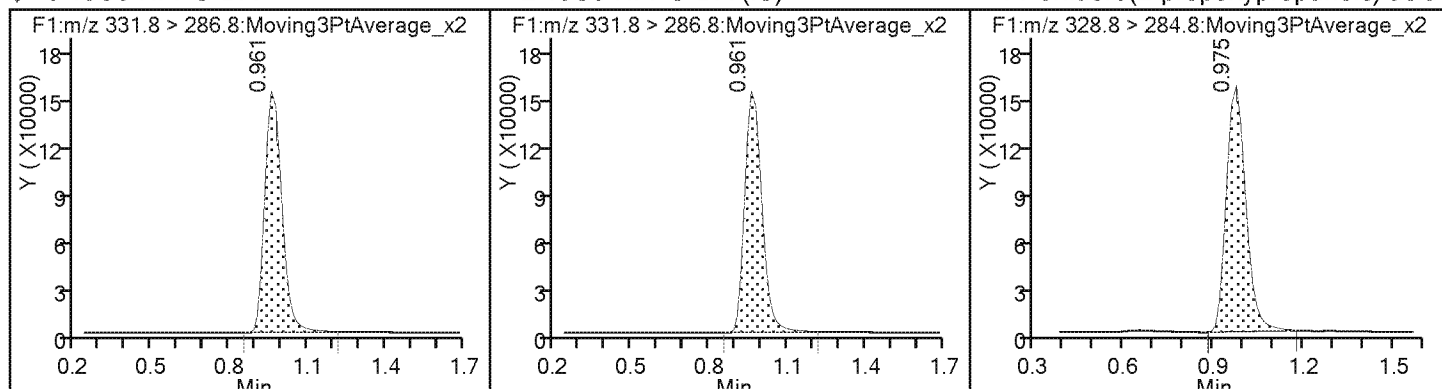

TestAmerica Denver  
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25008.d  
Lims ID: LCS 280-402648/2-A  
Client ID:  
Sample Type: LCS  
Inject. Date: 25-Jan-2018 10:44:48 ALS Bottle#: 12 Worklist Smp#: 5  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: LCS280-402648/2-A  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:16 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyera

Date: 25-Jan-2018 15:25:42

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	9.10	91.03



FORM I  
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: <u>TestAmerica Denver</u>	Job No.: <u>280-105698-1</u>
SDG No.: _____	
Client Sample ID: _____	Lab Sample ID: <u>LCSD 280-402648/3-A</u>
Matrix: <u>Water</u>	Lab File ID: <u>hfpo718A25009.d</u>
Analysis Method: <u>8321A</u>	Date Collected: _____
Extraction Method: <u>3535</u>	Date Extracted: <u>01/24/2018 15:20</u>
Sample wt/vol: <u>250 (mL)</u>	Date Analyzed: <u>01/25/2018 10:48</u>
Con. Extract Vol.: <u>5 (mL)</u>	Dilution Factor: <u>1</u>
Injection Volume: <u>20 (uL)</u>	GC Column: <u>Synergi Hydro</u> ID: _____
% Moisture: _____	GPC Cleanup: (Y/N) <u>N</u>
Analysis Batch No.: <u>402806</u>	Units: <u>ug/L</u>

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	0.190		0.010	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	90		50-200

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25009.d  
Lims ID: LCSD 280-402648/3-A  
Client ID:  
Sample Type: LCSD  
Inject. Date: 25-Jan-2018 10:48:03 ALS Bottle#: 13 Worklist Smp#: 6  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: LCSD280-402648/3-A  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:16 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyer

Date: 25-Jan-2018 15:25:45

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 &gt; 286.8 0.961 0.961 0.0 1.000 659150 9.01 1329

\* 2 13C3 HFPO-DA (IS)

331.8 &gt; 286.8 0.961 0.961 0.0 659150 10.0 1329

1 Perfluoro(2-propoxypropanoic) acid

328.8 &gt; 284.8 0.975 0.988 -0.013 1.000 648297 9.50 318

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25009.d

Injection Date: 25-Jan-2018 10:48:03

Instrument ID: LC\_LCMS7

Lims ID: LCSD 280-402648/3-A

Client ID:

Operator ID: JBH

ALS Bottle#: 13

Worklist Smp#: 6

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

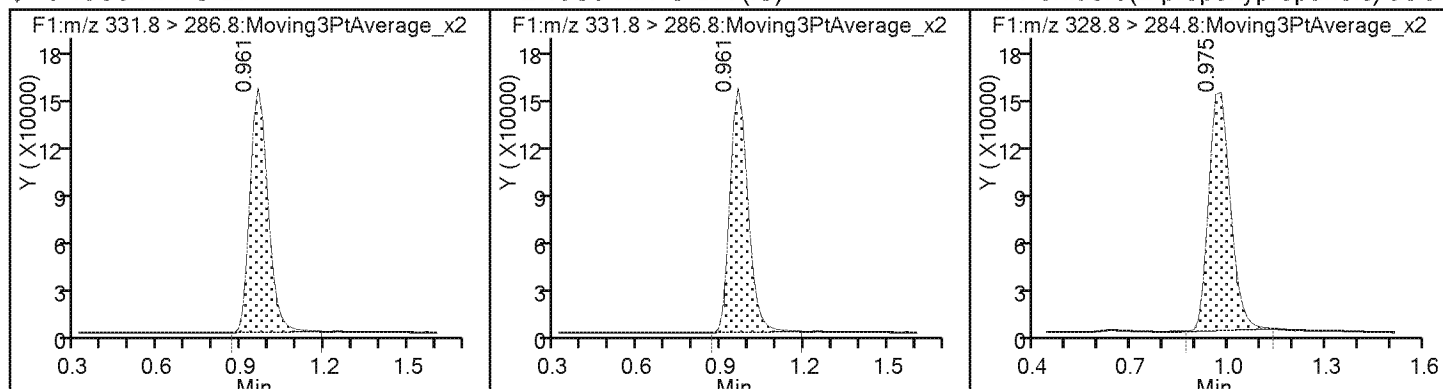

TestAmerica Denver  
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25009.d  
Lims ID: LCSD 280-402648/3-A  
Client ID:  
Sample Type: LCSD  
Inject. Date: 25-Jan-2018 10:48:03 ALS Bottle#: 13 Worklist Smp#: 6  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: LCSD280-402648/3-A  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:16 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyera

Date: 25-Jan-2018 15:25:45

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	9.01	90.12

FORM I  
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: <u>TestAmerica Denver</u>	Job No.: <u>280-105698-1</u>
SDG No.: _____	
Client Sample ID: _____	Lab Sample ID: <u>LLCS 280-402648/4-A</u>
Matrix: <u>Water</u>	Lab File ID: <u>hfpo718A25010.d</u>
Analysis Method: <u>8321A</u>	Date Collected: _____
Extraction Method: <u>3535</u>	Date Extracted: <u>01/24/2018 15:20</u>
Sample wt/vol: <u>250 (mL)</u>	Date Analyzed: <u>01/25/2018 10:51</u>
Con. Extract Vol.: <u>5 (mL)</u>	Dilution Factor: <u>1</u>
Injection Volume: <u>20 (uL)</u>	GC Column: <u>Synergi Hydro</u> ID: _____
% Moisture: _____	GPC Cleanup: (Y/N) <u>N</u>
Analysis Batch No.: <u>402806</u>	Units: <u>ug/L</u>

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	0.0173		0.010	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	94		50-200

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25010.d  
Lims ID: LLCS 280-402648/4-A  
Client ID:  
Sample Type: LLCS  
Inject. Date: 25-Jan-2018 10:51:18 ALS Bottle#: 14 Worklist Smp#: 7  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: LLCS280-402648/4-A  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:16 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyera

Date: 25-Jan-2018 15:25:48

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 &gt; 286.8 0.961 0.961 0.0 1.000 687452 9.40 1582

\* 2 13C3 HFPO-DA (IS)

331.8 &gt; 286.8 0.961 0.961 0.0 687452 10.0 1582

1 Perfluoro(2-propoxypropanoic) acid

328.8 &gt; 284.8 0.961 0.988 -0.027 1.000 75248 0.8656 34.8 M

## QC Flag Legend

Review Flags

M - Manually Integrated

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25010.d

Injection Date: 25-Jan-2018 10:51:18

Instrument ID: LC\_LCMS7

Lims ID: LLCS 280-402648/4-A

Client ID:

Operator ID: JBH

ALS Bottle#: 14

Worklist Smp#: 7

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid (M)

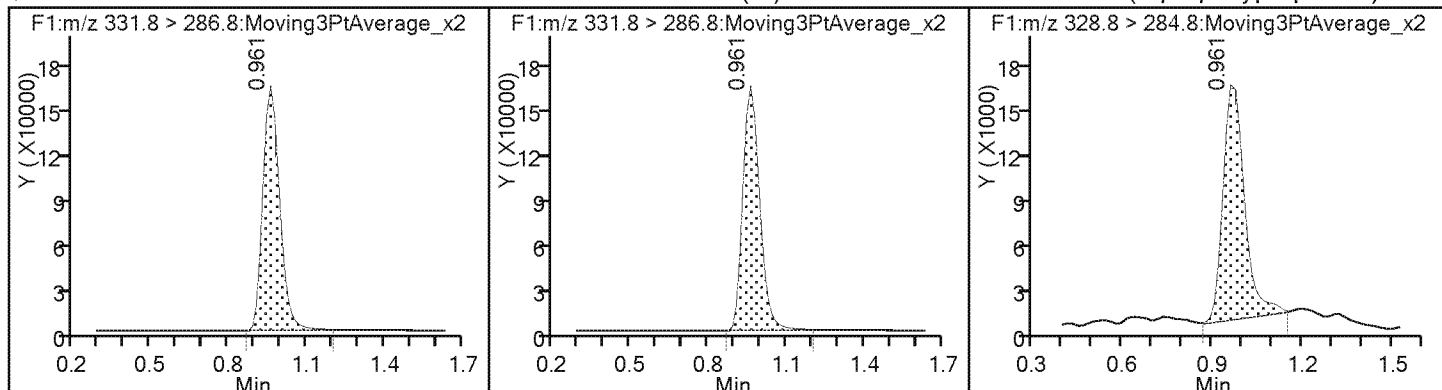

TestAmerica Denver  
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25010.d  
Lims ID: LLCS 280-402648/4-A  
Client ID:  
Sample Type: LLCS  
Inject. Date: 25-Jan-2018 10:51:18 ALS Bottle#: 14 Worklist Smp#: 7  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: LLCS280-402648/4-A  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:16 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyera

Date: 25-Jan-2018 15:25:48

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	9.40	93.99



## TestAmerica Denver


Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25010.d  
Injection Date: 25-Jan-2018 10:51:18 Instrument ID: LC\_LCMS7  
Lims ID: LLCs 280-402648/4-A  
Client ID:  
Operator ID: JBH ALS Bottle#: 14 Worklist Smp#: 7  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Method: HFPO Limit Group: LC - 8321A\_HFPO\_Du  
Column: Detector F1:MRM

## 1 Perfluoro(2-propoxypropanoic) acid, CAS: 13252-13-6

Signal: 1


RT: 0.96  
Area: 69234  
Amount: 0.779205  
Amount Units: ug/l

## Processing Integration Results

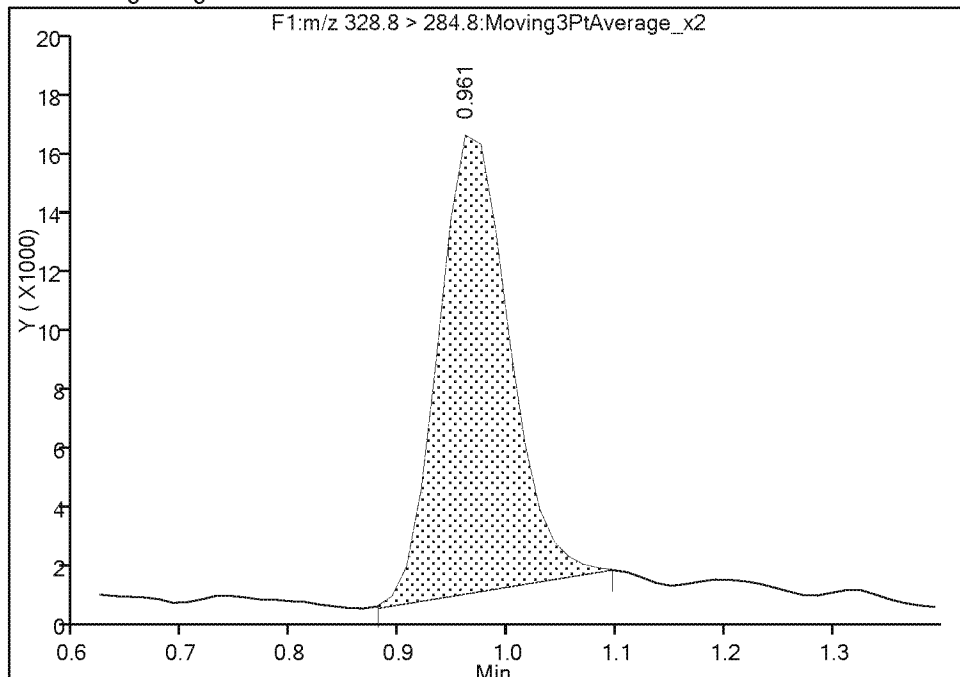

RT: 0.96  
Area: 75248  
Amount: 0.865642  
Amount Units: ug/l

## Manual Integration Results

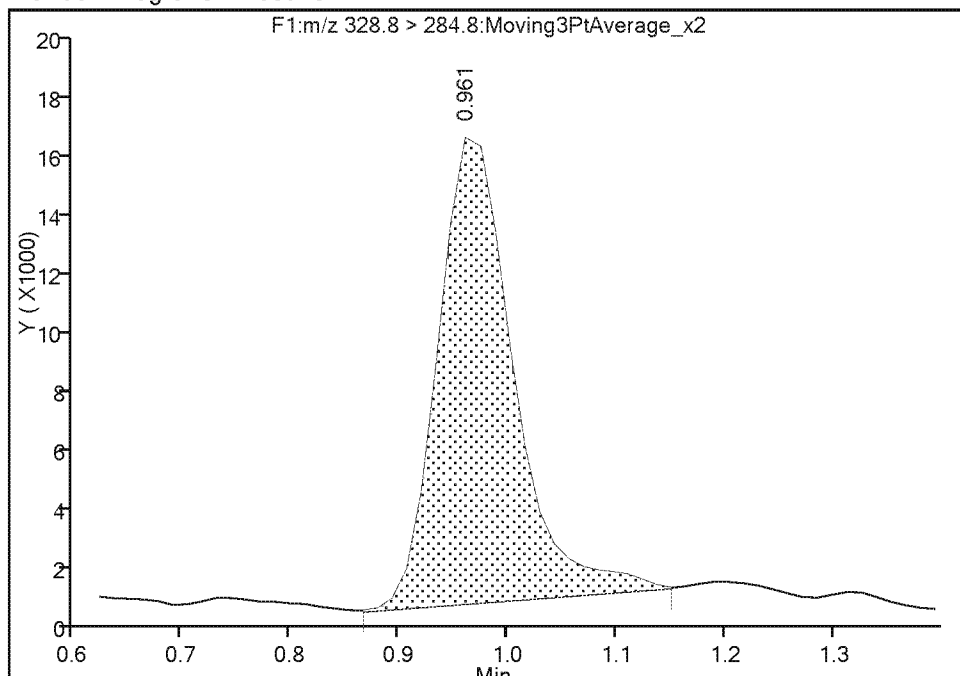

Reviewer: meyer, 25-Jan-2018 15:27:00

Audit Action: Manually Integrated

Audit Reason: Baseline

FORM I  
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: TestAmerica Denver Job No.: 280-105698-1  
 SDG No.: \_\_\_\_\_  
 Client Sample ID: \_\_\_\_\_ Lab Sample ID: DLCK 280-390728/12  
 Matrix: Water Lab File ID: hfpo717J10035.d  
 Analysis Method: 8321A Date Collected: \_\_\_\_\_  
 Extraction Method: \_\_\_\_\_ Date Extracted: \_\_\_\_\_  
 Sample wt/vol: 1(mL) Date Analyzed: 10/10/2017 10:04  
 Con. Extract Vol.: \_\_\_\_\_ Dilution Factor: 1  
 Injection Volume: 20(uL) GC Column: Synergi Hydro ID: \_\_\_\_\_  
 % Moisture: \_\_\_\_\_ GPC Cleanup: (Y/N) N  
 Analysis Batch No.: 390728 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	<0.50		0.50	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	102		50-200

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10035.d  
Lims ID: DLCK  
Client ID:  
Sample Type: DLCK  
Inject. Date: 10-Oct-2017 10:04:34 ALS Bottle#: 2 Worklist Smp#: 12  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: DLCK  
Misc. Info.: HFPO17J10  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 10-Oct-2017 12:51:51 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK005

First Level Reviewer: meyer

Date: 10-Oct-2017 11:51:31

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 &gt; 286.8 0.880 0.880 0.0 1.000 749614 10.2 480

\* 2 13C3 HFPO-DA (IS)

331.8 &gt; 286.8 0.880 0.880 0.0 749614 10.0 480

1 Perfluoro(2-propoxypropanoic) acid

328.8 &gt; 284.8 0.893 0.885 0.008 1.000 31104 0.1941 16.6

## Reagents:

HFPO\_CAL-1\_00031

Amount Added: 1.00

Units: mL

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10035.d

Injection Date: 10-Oct-2017 10:04:34

Instrument ID: LC\_LCMS7

Lims ID: DLCK

Client ID:

Operator ID: JBH

ALS Bottle#:

2

Worklist Smp#:

12

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

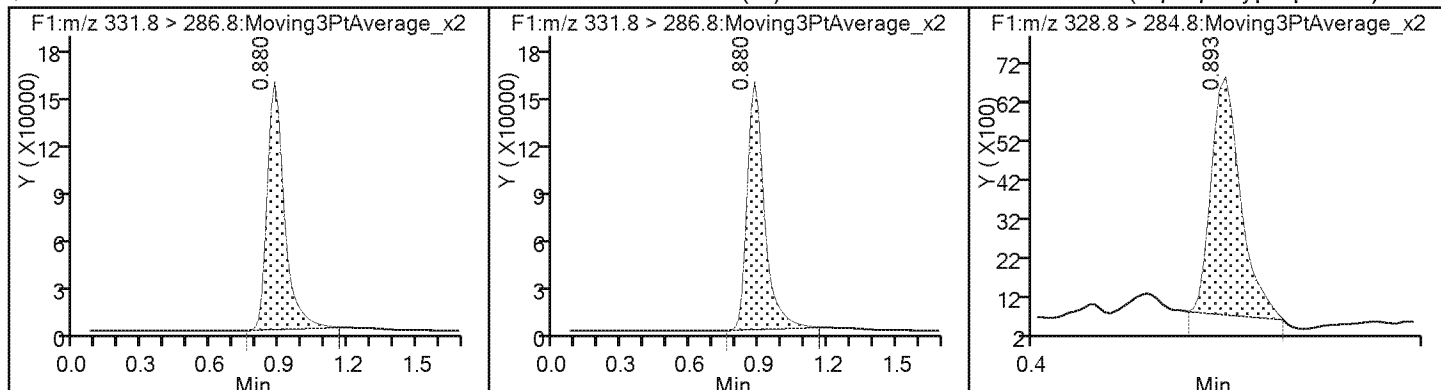

TestAmerica Denver  
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10035.d  
Lims ID: DLCK  
Client ID:  
Sample Type: DLCK  
Inject. Date: 10-Oct-2017 10:04:34 ALS Bottle#: 2 Worklist Smp#: 12  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: DLCK  
Misc. Info.: HFPO17J10  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 10-Oct-2017 12:51:51 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK005

First Level Reviewer: meyera

Date: 10-Oct-2017 11:51:31

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	10.2	102.48

Lab Name: <u>TestAmerica Denver</u>	Job No.: <u>280-105698-1</u>
SDG No.: _____	
Client Sample ID: <u>FAY-D-6394CHKFT-W1-1-0122</u> <u>18 MS</u>	Lab Sample ID: <u>280-105698-1 MS</u>
Matrix: <u>Water</u>	Lab File ID: <u>hfpo718A25013.d</u>
Analysis Method: <u>8321A</u>	Date Collected: <u>01/22/2018</u> <u>16:19</u>
Extraction Method: <u>3535</u>	Date Extracted: <u>01/24/2018</u> <u>15:20</u>
Sample wt/vol: <u>252.1(mL)</u>	Date Analyzed: <u>01/25/2018</u> <u>11:01</u>
Con. Extract Vol.: <u>5(mL)</u>	Dilution Factor: <u>1</u>
Injection Volume: <u>20(uL)</u>	GC Column: <u>Synergi Hydro</u> ID: _____
% Moisture: _____	GPC Cleanup: (Y/N) <u>N</u>
Analysis Batch No.: <u>402806</u>	Units: <u>ug/L</u>

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	0.220		0.010	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	86		50-200

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25013.d  
Lims ID: 280-105698-B-1-A MS  
Client ID: FAY-D-6394CHKFT-W1-1-012218  
Sample Type: MS  
Inject. Date: 25-Jan-2018 11:01:04 ALS Bottle#: 17 Worklist Smp#: 10  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: 280-105698-B-1-AMS  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:16 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyer

Date: 25-Jan-2018 15:26:36

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 &gt; 286.8 0.961 0.961 0.0 1.000 625788 8.56 1449

\* 2 13C3 HFPO-DA (IS)

331.8 &gt; 286.8 0.961 0.961 0.0 625788 10.0 1449

1 Perfluoro(2-propoxypropanoic) acid

328.8 &gt; 284.8 0.961 0.988 -0.027 1.000 716909 11.1 212

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfp0718A25013.d

Injection Date: 25-Jan-2018 11:01:04

Instrument ID: LC\_LCMS7

Lims ID: 280-105698-B-1-A MS

Client ID: FAY-D-6394CHKFT-W1-1-012218

Operator ID: JBH

ALS Bottle#: 17

Worklist Smp#: 10

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

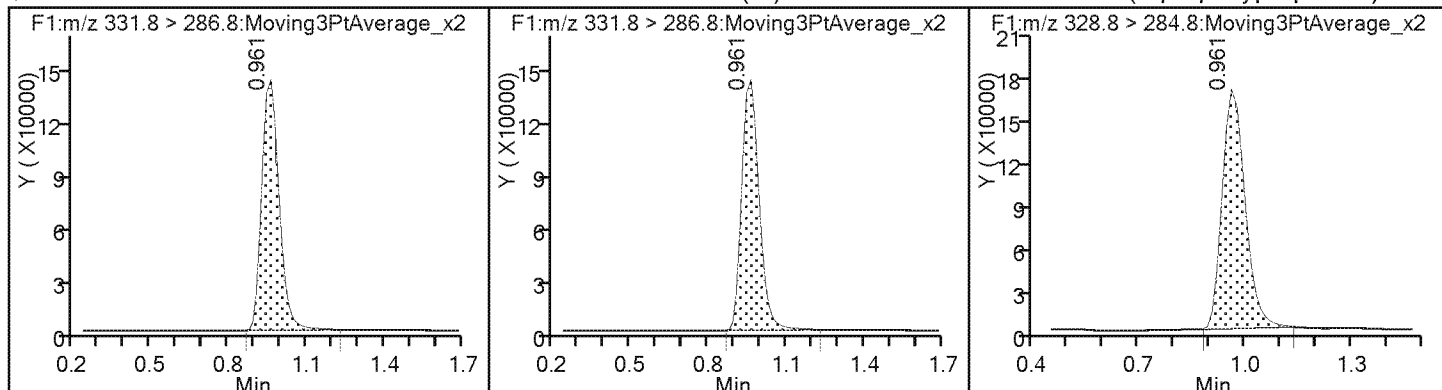



TestAmerica Denver  
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25013.d  
Lims ID: 280-105698-B-1-A MS  
Client ID: FAY-D-6394CHKFT-W1-1-012218  
Sample Type: MS  
Inject. Date: 25-Jan-2018 11:01:04 ALS Bottle#: 17 Worklist Smp#: 10  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: 280-105698-B-1-AMS  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:16 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyera

Date: 25-Jan-2018 15:26:36

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	8.56	85.55

## FORM I

Lab Name: TestAmerica Denver                      Job No.: 280-105698-1

SDG No.:

Client Sample ID: FAY-D-6394CHKFT-W1-1-0122    Lab Sample ID: 280-105698-1 DU  
18 DU

Matrix: Water Lab File ID: hfpo718A25012.d

Analysis Method: 8321A Date Collected: 01/22/2018 16:19

Extraction Method: 3535 Date Extracted: 01/24/2018 15:20

Sample wt/vol: 261.5 (mL) Date Analyzed: 01/25/2018 10:57

Con. Extract Vol.: 5(mL)                      Dilution Factor: 1

Injection Volume: 20(uL) GC Column: Synergi Hydro ID:

% Moisture: GPC Cleanup: (Y/N) N

Analysis Batch No.: 402806 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	0.0331		0.010	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	85		50-200

TestAmerica Denver  
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25012.d  
Lims ID: 280-105698-C-1-A DU  
Client ID: FAY-D-6394CHKFT-W1-1-012218  
Sample Type: DU  
Inject. Date: 25-Jan-2018 10:57:49 ALS Bottle#: 16 Worklist Smp#: 9  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: 280-105698-C-1-ADU  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:16 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyer

Date: 25-Jan-2018 15:26:18

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	-----------	-----------	-----------	----------	----------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 &gt; 286.8 0.961 0.961 0.0 1.000 619647 8.47 1456

\* 2 13C3 HFPO-DA (IS)

331.8 &gt; 286.8 0.961 0.961 0.0 619647 10.0 1456

1 Perfluoro(2-propoxypropanoic) acid

328.8 &gt; 284.8 0.961 0.988 -0.027 1.000 122268 1.73 36.7

## TestAmerica Denver

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25012.d

Injection Date: 25-Jan-2018 10:57:49

Instrument ID: LC\_LCMS7

Lims ID: 280-105698-C-1-A DU

Client ID: FAY-D-6394CHKFT-W1-1-012218

Operator ID: JBH

ALS Bottle#: 16

Worklist Smp#: 9

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A\_HFPO\_Du

\$ 3 13C3 HFPO-DA

\* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

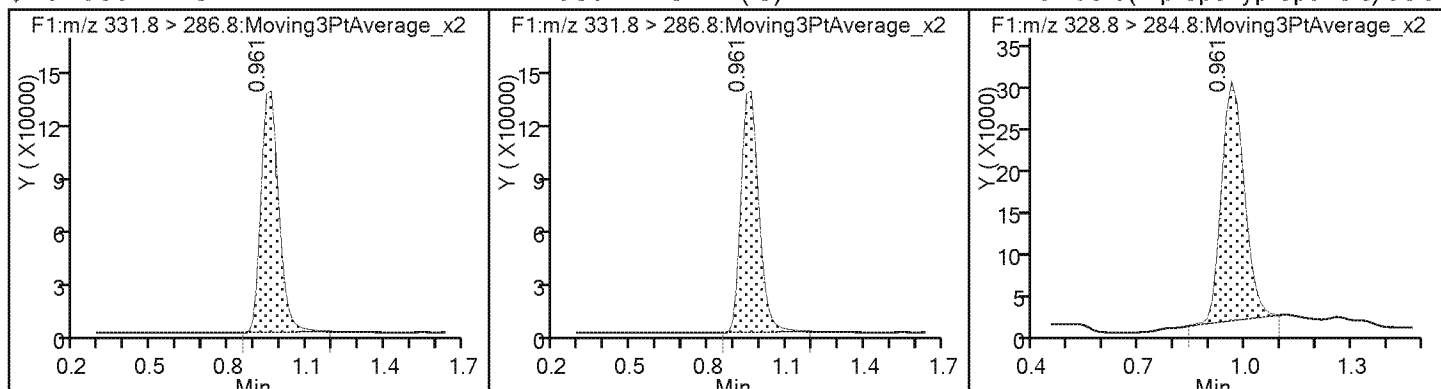

TestAmerica Denver  
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\hfpo718A25012.d  
Lims ID: 280-105698-C-1-A DU  
Client ID: FAY-D-6394CHKFT-W1-1-012218  
Sample Type: DU  
Inject. Date: 25-Jan-2018 10:57:49 ALS Bottle#: 16 Worklist Smp#: 9  
Injection Vol: 20.0 ul Dil. Factor: 1.0000  
Sample Info: 280-105698-C-1-ADU  
Misc. Info.: HFPO18A25  
Operator ID: JBH Instrument ID: LC\_LCMS7  
Method: \\ChromNA\Denver\ChromData\LC\_LCMS7\20180125-66715.b\HFPO.m  
Limit Group: LC - 8321A\_HFPO\_Du  
Last Update: 25-Jan-2018 15:31:16 Calib Date: 10-Oct-2017 09:58:07  
Integrator: Picker  
Quant Method: Internal/External Standard Quant By: Initial Calibration  
Last ICal File: \\ChromNA\Denver\ChromData\LC\_LCMS7\20171010-63483.b\hfpo717J10033.d  
Column 1 : Det: F1:MRM  
Process Host: XAWRK027

First Level Reviewer: meyera

Date: 25-Jan-2018 15:26:18

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	8.47	84.72

## LCMS ANALYSIS RUN LOG

Lab Name: TestAmerica DenverJob No.: 280-105698-1

SDG No.: \_\_\_\_\_

Instrument ID: LC\_LCMS7Start Date: 09/14/2017 14:40Analysis Batch Number: 387775End Date: 09/14/2017 16:05

LAB SAMPLE ID	CLIENT SAMPLE ID	DATE ANALYZED	DILUTION FACTOR	LAB FILE ID	COLUMN ID
STD001 280-387775/3 IC		09/14/2017 14:40	1	hfpo717I14052.d	Synergi Hydro
STD002 280-387775/4 IC		09/14/2017 14:43	1	hfpo717I14053.d	Synergi Hydro
STD003 280-387775/5 IC		09/14/2017 14:46	1	hfpo717I14054.d	Synergi Hydro
STD004 280-387775/6 IC		09/14/2017 14:49	1	hfpo717I14055.d	Synergi Hydro
STD005 280-387775/7 IC		09/14/2017 14:52	1	hfpo717I14056.d	Synergi Hydro
STD006 280-387775/8 IC		09/14/2017 14:55	1	hfpo717I14057.d	Synergi Hydro
STD007 280-387775/9 IC		09/14/2017 14:58	1	hfpo717I14058.d	Synergi Hydro
STD008 280-387775/10 IC		09/14/2017 15:01	1	hfpo717I14059.d	Synergi Hydro
ICB 280-387775/11		09/14/2017 15:04	1		Synergi Hydro
ZZZZZ		09/14/2017 15:07	1		Synergi Hydro
ICV 280-387775/13		09/14/2017 15:10	1	hfpo717I14062.d	Synergi Hydro
ZZZZZ		09/14/2017 15:13	1		Synergi Hydro
ZZZZZ		09/14/2017 15:16	1		Synergi Hydro
ZZZZZ		09/14/2017 15:19	1		Synergi Hydro
ZZZZZ		09/14/2017 15:22	1		Synergi Hydro
ZZZZZ		09/14/2017 15:25	2		Synergi Hydro
ZZZZZ		09/14/2017 15:28	1		Synergi Hydro
ZZZZZ		09/14/2017 15:31	1		Synergi Hydro
ZZZZZ		09/14/2017 15:34	1		Synergi Hydro
ZZZZZ		09/14/2017 15:38	1		Synergi Hydro
ZZZZZ		09/14/2017 15:41	1		Synergi Hydro
CCV 280-387775/24		09/14/2017 15:44	1		Synergi Hydro
ZZZZZ		09/14/2017 15:47	2		Synergi Hydro
ZZZZZ		09/14/2017 15:50	1		Synergi Hydro
ZZZZZ		09/14/2017 15:53	1		Synergi Hydro
ZZZZZ		09/14/2017 15:56	1		Synergi Hydro
ZZZZZ		09/14/2017 15:59	1		Synergi Hydro
ZZZZZ		09/14/2017 16:02	1		Synergi Hydro
CCV 280-387775/31		09/14/2017 16:05	1		Synergi Hydro

## LCMS ANALYSIS RUN LOG

Lab Name: TestAmerica DenverJob No.: 280-105698-1

SDG No.: \_\_\_\_\_

Instrument ID: LC\_LCMS7Start Date: 10/10/2017 09:35Analysis Batch Number: 390728End Date: 10/10/2017 11:19

LAB SAMPLE ID	CLIENT SAMPLE ID	DATE ANALYZED	DILUTION FACTOR	LAB FILE ID	COLUMN ID
STD001 280-390728/3 IC		10/10/2017 09:35	1	hfpo717J10026.d	Synergi Hydro
STD002 280-390728/4 IC		10/10/2017 09:38	1	hfpo717J10027.d	Synergi Hydro
STD003 280-390728/5 IC		10/10/2017 09:41	1	hfpo717J10028.d	Synergi Hydro
STD004 280-390728/6 IC		10/10/2017 09:45	1	hfpo717J10029.d	Synergi Hydro
STD005 280-390728/7 IC		10/10/2017 09:48	1	hfpo717J10030.d	Synergi Hydro
STD006 280-390728/8 IC		10/10/2017 09:51	1	hfpo717J10031.d	Synergi Hydro
STD007 280-390728/9 IC		10/10/2017 09:54	1	hfpo717J10032.d	Synergi Hydro
STD008 280-390728/10 IC		10/10/2017 09:58	1	hfpo717J10033.d	Synergi Hydro
ICB 280-390728/11		10/10/2017 10:01	1	hfpo717J10034.d	Synergi Hydro
DLCK 280-390728/12		10/10/2017 10:04	1	hfpo717J10035.d	Synergi Hydro
ICV 280-390728/13		10/10/2017 10:07	1	hfpo717J10036.d	Synergi Hydro
ZZZZZ		10/10/2017 10:11	1		Synergi Hydro
ZZZZZ		10/10/2017 10:14	1		Synergi Hydro
ZZZZZ		10/10/2017 10:17	1		Synergi Hydro
ZZZZZ		10/10/2017 10:20	1		Synergi Hydro
ZZZZZ		10/10/2017 10:23	1		Synergi Hydro
ZZZZZ		10/10/2017 10:27	1		Synergi Hydro
ZZZZZ		10/10/2017 10:30	1		Synergi Hydro
ZZZZZ		10/10/2017 10:33	1		Synergi Hydro
ZZZZZ		10/10/2017 10:36	1		Synergi Hydro
ZZZZZ		10/10/2017 10:40	1		Synergi Hydro
CCV 280-390728/24		10/10/2017 10:43	1	hfpo717J10047.d	Synergi Hydro
ZZZZZ		10/10/2017 10:46	1		Synergi Hydro
ZZZZZ		10/10/2017 10:49	1		Synergi Hydro
ZZZZZ		10/10/2017 10:53	1		Synergi Hydro
ZZZZZ		10/10/2017 10:56	1		Synergi Hydro
ZZZZZ		10/10/2017 10:59	1		Synergi Hydro
ZZZZZ		10/10/2017 11:02	1		Synergi Hydro
ZZZZZ		10/10/2017 11:06	1		Synergi Hydro
ZZZZZ		10/10/2017 11:09	1		Synergi Hydro
ZZZZZ		10/10/2017 11:12	1		Synergi Hydro
ZZZZZ		10/10/2017 11:16	1		Synergi Hydro
CCV 280-390728/35		10/10/2017 11:19	1		Synergi Hydro

## LCMS ANALYSIS RUN LOG

Lab Name: TestAmerica DenverJob No.: 280-105698-1

SDG No.: \_\_\_\_\_

Instrument ID: LC\_LCMS7Start Date: 01/25/2018 10:38Analysis Batch Number: 402806End Date: 01/25/2018 11:59

LAB SAMPLE ID	CLIENT SAMPLE ID	DATE ANALYZED	DILUTION FACTOR	LAB FILE ID	COLUMN ID
CCV 280-402806/3		01/25/2018 10:38	1	hfpo718A25006.d	Synergi Hydro
MB 280-402648/1-A		01/25/2018 10:41	1	hfpo718A25007.d	Synergi Hydro
LCS 280-402648/2-A		01/25/2018 10:44	1	hfpo718A25008.d	Synergi Hydro
LCSD 280-402648/3-A		01/25/2018 10:48	1	hfpo718A25009.d	Synergi Hydro
LLCS 280-402648/4-A		01/25/2018 10:51	1	hfpo718A25010.d	Synergi Hydro
280-105698-1		01/25/2018 10:54	1	hfpo718A25011.d	Synergi Hydro
280-105698-1 DU		01/25/2018 10:57	1	hfpo718A25012.d	Synergi Hydro
280-105698-1 MS		01/25/2018 11:01	1	hfpo718A25013.d	Synergi Hydro
CCV 280-402806/11		01/25/2018 11:04	1	hfpo718A25014.d	Synergi Hydro
280-105698-2		01/25/2018 11:07	1	hfpo718A25015.d	Synergi Hydro
280-105698-3		01/25/2018 11:10	1	hfpo718A25016.d	Synergi Hydro
280-105698-4		01/25/2018 11:14	1	hfpo718A25017.d	Synergi Hydro
280-105698-5		01/25/2018 11:17	1	hfpo718A25018.d	Synergi Hydro
280-105698-6		01/25/2018 11:20	1	hfpo718A25019.d	Synergi Hydro
280-105698-7		01/25/2018 11:23	1	hfpo718A25020.d	Synergi Hydro
280-105698-8		01/25/2018 11:27	1	hfpo718A25021.d	Synergi Hydro
CCV 280-402806/19		01/25/2018 11:30	1	hfpo718A25022.d	Synergi Hydro
280-105698-9		01/25/2018 11:33	1	hfpo718A25023.d	Synergi Hydro
280-105698-10		01/25/2018 11:36	1	hfpo718A25024.d	Synergi Hydro
ZZZZZ		01/25/2018 11:40	1		Synergi Hydro
ZZZZZ		01/25/2018 11:43	1		Synergi Hydro
ZZZZZ		01/25/2018 11:46	1		Synergi Hydro
ZZZZZ		01/25/2018 11:49	1		Synergi Hydro
ZZZZZ		01/25/2018 11:53	1		Synergi Hydro
ZZZZZ		01/25/2018 11:56	1		Synergi Hydro
CCV 280-402806/28		01/25/2018 11:59	1	hfpo718A25031.d	Synergi Hydro



## LCMS BATCH WORKSHEET

Lab Name: TestAmerica DenverJob No.: 280-105698-1

SDG No.: \_\_\_\_\_

Batch Number: 402648Batch Start Date: 01/24/18 15:20Batch Analyst: Mueller, Stacey KBatch Method: 3535Batch End Date: 01/24/18 20:35

Lab Sample ID	Client Sample ID	Method Chain	Basis	GrossWeight	TareWeight	InitialAmount	FinalAmount	HFPO I.S. 00007	HFPO Spike 00004
MB 280-402648/1		3535, 8321A				250 mL	5 mL	0.1 mL	
LCS 280-402648/2		3535, 8321A				250 mL	5 mL	0.1 mL	0.1 mL
LCSD 280-402648/3		3535, 8321A				250 mL	5 mL	0.1 mL	0.1 mL
LLCS 280-402648/4		3535, 8321A				250 mL	5 mL	0.1 mL	0.01 mL
280-105698-D-1	FAY-D-6394CHKFT-W1-1-012218	3535, 8321A	T	290.8 g	27.0 g	263.8 mL	5 mL	0.1 mL	
280-105698-C-1 DU	FAY-D-6394CHKFT-W1-1-012218	3535, 8321A	T	289.1 g	27.6 g	261.5 mL	5 mL	0.1 mL	
280-105698-B-1 MS	FAY-D-6394CHKFT-W1-1-012218	3535, 8321A	T	280.3 g	28.2 g	252.1 mL	5 mL	0.1 mL	0.1 mL
280-105698-B-2	FAY-D-6394CHKFT-W1-1-012218-D	3535, 8321A	T	278.3 g	27.6 g	250.7 mL	5 mL	0.1 mL	
280-105698-D-3	FAY-D-6246CHKFT-W1-1-012218	3535, 8321A	T	280.9 g	26.8 g	254.1 mL	5 mL	0.1 mL	
280-105698-A-4	FAY-D-318BOONE-W1-1-012218	3535, 8321A	T	272.6 g	27.5 g	245.1 mL	5 mL	0.1 mL	
280-105698-D-5	FAY-D-41BOONE-W1-1-012218	3535, 8321A	T	280.8 g	28.1 g	252.7 mL	5 mL	0.1 mL	
280-105698-D-6	FAY-D-FB-012218	3535, 8321A	T	288.3 g	27.9 g	260.4 mL	5 mL	0.1 mL	
280-105698-C-7	FAY-D-7145BUTLE-W1-1-012218	3535, 8321A	T	283.5 g	26.4 g	257.1 mL	5 mL	0.1 mL	
280-105698-B-8	FAY-D-1515SCLLY-W1-1-012218	3535, 8321A	T	281.3 g	27.2 g	254.1 mL	5 mL	0.1 mL	
280-105698-A-9	FAY-D-7396SALIE-W1-1-012218	3535, 8321A	T	287.6 g	29.0 g	258.6 mL	5 mL	0.1 mL	
280-105698-A-10	FAY-D-7012NC87H-W1-1-012218	3535, 8321A	T	275.0 g	27.3 g	247.7 mL	5 mL	0.1 mL	

The pound sign (#) in the amount added field denotes that the reagent was used undiluted. All calculations are performed using the stated concentration for this reagent.

8321A

Page 1 of 2

## LCMS BATCH WORKSHEET

Lab Name: TestAmerica DenverJob No.: 280-105698-1

SDG No.: \_\_\_\_\_

Batch Number: 402648Batch Start Date: 01/24/18 15:20Batch Analyst: Mueller, Stacey KBatch Method: 3535Batch End Date: 01/24/18 20:35

Batch Notes	
Acid ID	2% Formic Aci_00139
Acid Name	2% Formic Acid
Balance ID	24350888
Batch Comment	Reviewer: AMB
First End time	1700
H2O ID	HPLC_Water_00845/846
Pipette ID	m2, SPE-1 + syringe (LLCS)
Reagent ID	10% NH4OH
Reagent Lot Number	10% NH4OH_00116
Solvent Lot #	Methanol_00188
Solvent Name	Methanol
SOP Number	DV-OP-0019
SPE Cartridge Type	STRATA-X-AW (8B S038 FCH)
Solid Phase Extraction Disk ID	S308-0077
First Start time	1540

Basis	Basis Description
T	Total/NA

The pound sign (#) in the amount added field denotes that the reagent was used undiluted. All calculations are performed using the stated concentration for this reagent.

8321A

Page 2 of 2



**Reagent ID: HFPO\_CAL-6\_00078**

Description: level8  
No. of Bottles: 1  
Storage Location: LCMS  
Reagent Volume: 1.000 mL  
Creation Date: 01/12/2018  
Open Date:  
Container(s): 4921808  
Comment: level-8

Expiration Date: 01/28/2018  
Laboratory: TestAmerica Denver  
Prepared By: Meyer, Andrew GC  
Solvent: 80:20 Methanol : H2O  
Solvent Lot: 00016

### Reagent Analyte Information

Analyte	Source ID	Source Exp. Date	Source Conc.	Source Conc. Units	Final Conc.	Final Conc. Units
13C3 HFPO-DA	HFPO I.S._00004	08/28/2018	0.50000	ug/mL	10.00000	ug/L
13C3 HFPO-DA (IS)	HFPO I.S._00004	08/28/2018	0.50000	ug/mL	10.00000	ug/L
Perfluoro(2-propoxypropenol) acid	HFPO Spike_00004	10/30/2018	0.50000	ug/mL	10.00000	ug/L

### Source Reagents

Reagent	Description	Type	Expiration	Vendor	Vendor Lot #	Vendor Cat Lot #	Volume Used	Volume Units
HFPO I.S._00004	Internal Standard for HFPO 0.5ug/ml		08/28/18				20.00000	uL
HFPO Spike_00004	HFPO LCS/Calibration Spike 0.5ug/ml		10/30/18				20.00000	uL

*Andrew Meyer*  
01/26/2018



**Reagent ID: HFPO\_CAL-5\_00078**

Description: level5  
No. of Bottles: 1  
Storage Location: LCMS  
Reagent Volume: 1.000 mL  
Creation Date: 01/12/2018  
Open Date:  
Container(s): 4821807  
Comment: level-5

Expiration Date: 01/28/2018  
Laboratory: TestAmerica Denver  
Prepared By: Mayer, Andrew GC  
Solvent: 80:20 Methanol : H2O  
Solvent Lot: 00016

### Reagent Analyte Information

Analyte	Source ID	Source Exp. Date	Source Conc.	Source Conc. Units	Final Conc.	Final Conc. Units
13CS HFPO-DA	HFPO I.S._00004	08/28/2018	0.50000	ug/mL	10.00000	ug/L
13CS HFPO-DA (18)	HFPO I.S._00004	08/28/2018	0.50000	ug/mL	10.00000	ug/L
Perfluoro(2-propoxypropanoic) acid	HFPO Spike_00004	10/30/2018	0.50000	ug/mL	5.00000	ug/L

### Source Reagents

Reagent	Description	Type	Expiration	Vendor	Vendor Lot #	Vendor Cat Lot #	Volume Used	Volume Units
HFPO I.S._00004	Internal Standard for HFPO 0.5ug/ml		08/28/18				20.00000	uL
HFPO Spike_00004	HFPO LCS/Calibration Spike 0.5ug/ml		10/30/18				10.00000	uL

*Andrew Mayer*  
*01/12/2018*

# Shipping and Receiving Documents

**Estimotech America Denver**  
1955 Yarrow Street  
Arvada, CO 80002  
Phone (303) 735-0100 Fax (303) 735-0101

105683 Denied COL received -93

**TestAmerica**

WJ 44485-45. 2013. *Journal of Cellular Biochemistry* 109: 2560-2568.

<b>Client Information</b> Client Contact: <u>Mr. Michael Aucoin</u> Company: <u>The Chemours Company FC, LLC</u> Address: <u>4051 Ogletown Road, Suite 300</u> City: <u>Newark</u> State: <u>DE</u> Zip: <u>19713</u> Phone: <u>302.781.5873</u> Email: <u>michael.aucoin@aeocom.com</u> Project Name: <u>FAY-2018 Residential Sampling</u> Site: <u></u>		Samples: <u>AM, NV</u> Phone: <u>734-600-5746</u> E-Mail: <u>Michelle.johnston@lestameric.com</u> Job # <u>2018-01</u>		Lab Pkg: <u>Johnston, Michelle</u> E-Mail: <u>Michelle.johnston@lestameric.com</u> Date: <u>01/22/18</u>	
Due Date Requested: <u>10 Business Days</u> PO #: <u>LEJO-67048/64201000-2231QS1000</u> Project #: <u>28016904</u> SSOR#: <u></u>		Analysis Requested: <u></u> Total Number of Containers: <u></u> Special Instructions/Note: <u>Hold all remaining volumes as Reducers</u>			
Sample Identification Sample Date: <u>10/21/18</u> Sample Time: <u>1621 G</u> Sample Type: <u>G</u> Matrix: <u>W</u> Sample Date: <u>10/21/18</u> Sample Time: <u>1621 G</u> Sample Type: <u>G</u> Matrix: <u>W</u>		Possible Hazard Identification: <input checked="" type="checkbox"/> Non-Hazard <input type="checkbox"/> Flammable <input type="checkbox"/> Skin Irritant <input type="checkbox"/> Poison B <input type="checkbox"/> Unknown <input type="checkbox"/> Radiological Deliverable Requested: 1, II, III, IV, Other (specify) Level IV			
Empty Kit Requisitioned by: <u></u> Requisitioned by: <u></u> Requisitioned by: <u></u>		Sample Disposal (A fee may be assessed if samples are retained longer than 1 month) <input type="checkbox"/> Return To Client <input type="checkbox"/> Disposal By Lab <input type="checkbox"/> Archive For Months: <u></u> Special Instructions/QC Requirements: <u></u>			
Date/Time: <u>01/22/18 / 1600</u> Date/Time: <u></u> Date/Time: <u></u>		Date/Time: <u>01/22/18 0930</u> Date/Time: <u></u> Date/Time: <u></u>			
Requested by: <u></u> Requisitioned by: <u></u> Requisitioned by: <u></u>		Received by: <u></u> Received by: <u></u> Received by: <u></u>			
Company: <u></u> Company: <u></u> Company: <u></u>		Company: <u></u> Company: <u></u> Company: <u></u>			
Custody Seals Intact: <u></u> A Yes A No		Custody Seal No.: <u></u>			

## Chain of Custody Record

<b>Client Information</b> Mr. Michael Aucoin Company: The Chemours Company FC, LLC Address: c/o AECOM 4051 Ogletown Road, Suite 300 City: Newark State: DE 19713 Phone: 302.781.5973 Email: michael.aucoin@aecom.com Project Name: FAY-2018 Residential Sampling Site:		Shipper: KS, CM Date Recd: Johnston, Michelle E-Mail: michelle.johnston@testamerica.com Carrier Tracking Ref: FedEx		COC No: 1 of 2 Page: of Job #:		
<b>Analysis Requested</b> Due Date Requested: 10 Business Days TAT Requested (days): 10 Business Days TO: LBO-87048/84201000-2231051000 WQ: R Project #: 28016004 SOW#:						
<b>Sample Identification</b>						
Sample Identification	Sample Date	Sample Time	Sample Type (IC=Comp, G=grab)	Matrix (Inorganic, Organic, Suspended, Dissolved, Other)	Preservation Code	Matrix
FAY-D-6394CHKFT-WI-1-012218	01-22-18	1619	G	W	W	W
FAY-D-6394CHKFT-WI-1-012218-D	01-22-18	1619	G	W	W	W
FAY-D-6394CHKFT-WI-1-012218-REP	01-22-18	1619	G	W	W	W
FAY-D-6394CHKFT-WI-1-012218-MS	01-22-18	1619	G	W	W	W
FAY-D-6394CHKFT-WI-1-012218	01-22-18	1456	G	W	W	W
FAY-D-318BOONE-WI-1-012218	01-22-18	1415	G	W	W	W
FAY-D-418BOONE-WI-1-012218	01-22-18	1349	G	W	W	W
FAY-D-FB-012218	01-22-18	0730	G	W	W	W
FAY-D-7145BVTLE-WI-1-012218	01-22-18	1403	G	W	W	W
FAY-D-1515SCLLY-WI-1-012218	01-22-18	1423	G	W	W	W
FAY-D-7396SALJE-WI-1-012218	01-22-18	1539	G	W	W	W
Possible Hazard Identification <input type="checkbox"/> Non-Hazard <input type="checkbox"/> Flammable <input type="checkbox"/> Skin Irritant <input type="checkbox"/> Poison B <input type="checkbox"/> Unknown <input type="checkbox"/> Radiological Deliverable Requested: I, II, III, IV, Other (specify) Level IV						
Sample Disposal (A fee may be assessed if samples are retained longer than 1 month) <input type="checkbox"/> Return To Client <input type="checkbox"/> Disposal By Lab <input type="checkbox"/> Archive For Months						
Special Instructions/OC Requirements:						
Special Instructions/Notes: Hold all Volumes as Retains						
Total Number of Containers: 4						
Preservation Codes: A - HCL B - NaOH C - Zn Acetate D - Nitric Acid E - Nitric Acid F - NaOH G - Ascorbic Acid H - Ascorbic Acid I - Bz J - GI Water K - EDTA L - EDA M - HCL N - None O - Ascorbic Acid P - NaOH Q - NaOH R - NaOH S - H2SO4 T - TSP Desiccant U - Ascorbic Acid V - MCAA W - pH 4.5 X - Other (specify)						
Received by: Taylor Paswith Date Recd: 01-22-18 1700 Company: Paswith						
Custody Seal No.: Yes No						

## Chain of Custody Record

[illegible]



## Login Sample Receipt Checklist

Client: Chemours Company FC, LLC The

Job Number: 280-105698-1

Login Number: 105698

List Source: TestAmerica Denver

List Number: 1

Creator: True, Joshua A

Question	Answer	Comment
Radioactivity wasn't checked or is $\leq$ background as measured by a survey meter.	N/A	
The cooler's custody seal, if present, is intact.	True	Refer to Job Narrative for details.
Sample custody seals, if present, are intact.	True	
The cooler or samples do not appear to have been compromised or tampered with.	True	
Samples were received on ice.	True	
Cooler Temperature is acceptable.	True	
Cooler Temperature is recorded.	True	
COC is present.	True	
COC is filled out in ink and legible.	True	
COC is filled out with all pertinent information.	True	
Is the Field Sampler's name present on COC?	N/A	
There are no discrepancies between the containers received and the COC.	False	Refer to Job Narrative for details.
Samples are received within Holding Time (excluding tests with immediate HTs)	True	
Sample containers have legible labels.	True	
Containers are not broken or leaking.	True	
Sample collection date/times are provided.	True	
Appropriate sample containers are used.	True	
Sample bottles are completely filled.	True	
Sample Preservation Verified.	N/A	
There is sufficient vol. for all requested analyses, incl. any requested MS/MSDs	True	
Containers requiring zero headspace have no headspace or bubble is $<6\text{mm}$ (1/4").	N/A	
Multiphasic samples are not present.	True	
Samples do not require splitting or compositing.	True	
Residual Chlorine Checked.	N/A	