

ANALYTICAL REPORT

Job Number: 320-35355-1

Job Description: FAY-2018 Residential Sampling

For:

Chemours Company FC, LLC
c/o AECOM
Sabre Building, Suite 300
4051 Ogletown Road
Newark, DE 19713

Attention: Michael Aucoin



Approved for release.
Michelle A Johnston
Project Manager II
3/29/2018 1:05 PM

Michelle A Johnston, Project Manager II
4955 Yarrow Street, Arvada, CO, 80002
(303)736-0110
michelle.johnston@testamericainc.com
03/29/2018
Revision: 1

cc: Barbara McGraw
Kelly Rinehimer

The test results in this report relate only to the samples in this report and meet all requirements of NELAC, with any exceptions noted. Pursuant to NELAP, this report shall not be reproduced except in full, without the written approval of the laboratory. All questions regarding this report should be directed to the TestAmerica Denver Project Manager.

The Lab Certification ID# is 4025.

Reporting limits are adjusted for sample size used, dilutions and moisture content if applicable.

TestAmerica Laboratories, Inc.

TestAmerica Sacramento 880 Riverside Parkway, West Sacramento, CA 95605
Tel (916) 373-5600 Fax (916) 372-1059 www.testamericainc.com

Table of Contents

Cover Title Page	1
Data Summaries	4
Definitions	4
Case Narrative	5
Detection Summary	7
Client Sample Results	19
Default Detection Limits	21
Surrogate Summary	22
QC Sample Results	23
QC Association	25
Chronicle	26
Certification Summary	28
Method Summary	29
Sample Summary	30
Manual Integration Summary	31
Reagent Traceability	33
COAs	39
Organic Sample Data	65
LCMS	65
8321A_HFPO_Du	65
8321A_HFPO_Du QC Summary	66
8321A_HFPO_Du Sample Data	73
Standards Data	116
8321A_HFPO_Du ICAL Data	116
8321A_HFPO_Du CCAL Data	159
Raw QC Data	177

Table of Contents

8321A_HFPO_Du Tune Data	177
8321A_HFPO_Du Blank Data	182
8321A_HFPO_Du LCS/LCSD Data	190
8321A_HFPO_Du MS/MSD Data	207
8321A_HFPO_Du Duplicate/Tripleate Data	211
8321A_HFPO_Du Run Logs	216
8321A_HFPO_Du Prep Data	219
Shipping and Receiving Documents	223
Client Chain of Custody	224
Sample Receipt Checklist	228

Definitions/Glossary

Client: Chemours Company FC, LLC The
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 320-35355-1

Glossary

Abbreviation	These commonly used abbreviations may or may not be present in this report.
%	Listed under the "D" column to designate that the result is reported on a dry weight basis
%R	Percent Recovery
CFL	Contains Free Liquid
CNF	Contains No Free Liquid
DER	Duplicate Error Ratio (normalized absolute difference)
Dil Fac	Dilution Factor
DL	Detection Limit (DoD/DOE)
DL, RA, RE, IN	Indicates a Dilution, Re-analysis, Re-extraction, or additional Initial metals/anion analysis of the sample
DLC	Decision Level Concentration (Radiochemistry)
EDL	Estimated Detection Limit (Dioxin)
LOD	Limit of Detection (DoD/DOE)
LOQ	Limit of Quantitation (DoD/DOE)
MDA	Minimum Detectable Activity (Radiochemistry)
MDC	Minimum Detectable Concentration (Radiochemistry)
MDL	Method Detection Limit
ML	Minimum Level (Dioxin)
NC	Not Calculated
ND	Not Detected at the reporting limit (or MDL or EDL if shown)
PQL	Practical Quantitation Limit
QC	Quality Control
RER	Relative Error Ratio (Radiochemistry)
RL	Reporting Limit or Requested Limit (Radiochemistry)
RPD	Relative Percent Difference, a measure of the relative difference between two points
TEF	Toxicity Equivalent Factor (Dioxin)
TEQ	Toxicity Equivalent Quotient (Dioxin)

CASE NARRATIVE

Client: The Chemours Company FC, LLC

Project: FAY-2018 Residential Sampling

Report Number: 320-35355-1

With the exceptions noted as flags or footnotes, standard analytical protocols were followed in the analysis of the samples and no problems were encountered or anomalies observed. In addition all laboratory quality control samples were within established control limits, with any exceptions noted below. Each sample was analyzed to achieve the lowest possible reporting limit within the constraints of the method. In some cases, due to interference or analytes present at high concentrations, samples were diluted. For diluted samples, the reporting limits are adjusted relative to the dilution required.

Calculations are performed before rounding to avoid round-off errors in calculated results.

All holding times were met and proper preservation noted for the methods performed on these samples, unless otherwise detailed in the individual sections below.

For samples requiring analysis at a dilution, the dilution factor has been multiplied by the Method Detection Limit (MDL) for each analyte and evaluated versus the project-specific reporting limit (PSRL). If the obtained value is below the PSRL, then the PSRL is preserved as the reporting limit for the diluted result, otherwise, the obtained value becomes the reporting limit. This is done in order to maintain the PSRL to meet project requirements at the request of the client and to report the lowest possible RL for each analyte.

Revision - 3/29/2018

The certificate form for sample FAY-D-661PINEA-W1-1-012418 (320-35355-4) was revised to include the trailing zero in the final result.

Receipt

The samples were received on 1/25/2018 9:30 AM; the samples arrived in good condition, properly preserved and, where required, on ice. The temperature of the cooler at receipt was 0.6°C.

Receipt Exceptions

The laboratory logged all sample IDs per the revised chain of custody submitted to the laboratory by the client via email on 1/26/2018. The revised and original chain of custody are included in the final report.

Due to client shipping error, the parent sample volume for FAY-D-7166BUTLE-W1-1-012418 (320-35355-9) was received at TestAmerica Denver on 1/25/2018 and was originally logged under job 280-105815. Per client instructions received on 1/26/2018, the sample volume was re-logged to be included with sample volume from the same sample location designated for laboratory MS/REP QC, originally received at TestAmerica Sacramento on 1/25/2018, and logged under this job, 320-35355.

The requested HFPO-DA analyses were subcontracted to TestAmerica's Denver laboratory.

No other anomalies were observed during sample receipt.

Standards

Analytical standards were prepared using the acid form of the compound Perfluoro(2-propoxypropanoic) acid (HFPO-DA).

The surrogate compound, 13C3 HFPO-DA was introduced at the extraction step and was used as an internal standard for quantitation of HFPO-DA. The concentration of the surrogate spike is 0.2ug/L in water samples or 50ug/kg in soil samples.

Sample Extraction and Analysis

The samples presented in this report were extracted for the target analyte by TestAmerica Denver's SOP DV-OP-0019, Rev. 8 and analyzed for the target analyte by TestAmerica Denver's SOP DV-LC-0012, Rev. 14, with the exceptions of the items indicated in the DuPont QAS. Sample FAY-D-7166BUTLE-W1-1-012418 (320-35355-9) was chosen to be analyzed as a duplicate and also to be spiked with the target analyte.

For water samples a 250mL aliquot of each sample is extracted using solid phase extraction technique with methanol conditioned Weak Anion Exchange cartridges. Each sample is spiked with the internal standard/surrogate, prior to extraction. After the sample is passed through the cartridge, the analytes are eluted with 2%Formic Acid, 6mLs of HPLC grade MeOH and then with 4mL of 10% ammonium hydroxide in methanol. The final volume is brought to 5mL using reagent water and the extract is analyzed by LC/MS/MS.

The target analyte is separated from other components on a high-performance liquid chromatography (HPLC) C18 column with a mobile phase mixture of water containing 0.1% ammonium acetate and methanol. The mass spectrometer detector is operated in the electrospray (ESI) negative ion mode. The instrument is calibrated at 7 concentration levels (0.2, 0.5, 1.0, 2.0, 5.0, 10 and 20ug/L). The target analyte is detected as the perfluoro(2-propoxypropanoic) acid with the parent ion of 328.8 amu. The daughter ions used for analysis by LC/MS/MS are at 284.8 amu. The ratio of the peak areas to the two ions must be $\pm 20\%$ of the ion ratios in the mid-point ICAL for qualitative identification. Sample results are quantitated using the internal standard dilution.

Tuning and Calibration

The instrument is tuned with a solution of the target analyte such that mass assignments are within ± 0.5 amu of the daughter ions. The

instrument is calibrated with seven concentration levels from 0.2ug/L to 20ug/L. Linear regression ($y=ax+b$) or quadratic functions ($y=ax^2+bx+c$) are used with a correlation coefficient or coefficient of determination ≥ 0.990 . Following initial calibration (ICAL), an initial calibration blank (ICB) is tested, which consists of methanol spiked with the surrogate. The result for the target analyte must be less than one half the reporting limit (RL) to proceed.

Next an initial calibration verification (ICV) standard is tested. This is a mid-level concentration standard from a different vendor from the ICAL standard. If a different vendor is not available then, a different lot number from the same vendor is used. The ICV must be within 80-120% of the true value.

The quantitation limit verification standard is a standard from the same source as the ICAL tested run at the RL level to determine accuracy near the detection limit. This recovery must be within 70-130%.

Continuing calibration verification (CCV) standards are tested every 10 injections and are from the same source as the ICAL and are at mid-level concentration. The recovery of the CCVs must be 70-130% or recalibration is necessary.

Method QC Samples

The Method Blank is processed reagent water spiked with internal standard and prepared with each batch of 20 samples of the same matrix. All samples in the batch are processed at the same time and with the same reagents. The method blank must be less than the LOD or associated batch samples must be re-extracted and reanalyzed.

Each batch is prepared with a low- and a mid-level concentration spike Laboratory Control Samples (LCS). The recoveries of these samples must be within 70-130% or associated batch samples must be re-extracted and reanalyzed. If the recovery is biased high and samples are non-detect, results can be reported without re-extraction.

Calculations

Sample Result Calculation

For internal standard quantitation,

$$\text{HFPO-DA Response} = \text{Area of HFPO-DA} * 13C3 \text{ HFPO-DA concentration / area of } 13C3 \text{ HFPO-DA}$$

$$\text{Concentration in waters, ug/L} = (\text{Cex Vt})/(\text{Vo})$$

Where:

Cex = Concentration measured in sample extract from the target analyte response (ng/mL)

Vt = Volume of total extract (mL)

Vo = Volume of water extracted (mL)

2. Percent Recovery Calculation

$$\text{Spike Recovery} = (\text{SSR}-\text{SR})/(\text{SA})x100\%$$

Where:

SSR = Spike sample result

SR = Sample result

SA = Spike added

3. Relative Percent Difference Calculation

$$\text{RPD} = (\text{SR} - \text{DR})/(1/2(\text{SR}+\text{DR}))x100$$

Where:

SR = Sample result

DR = Duplicate result

HFPO-DA Analysis Anomalies

Samples FAY-D-7154BUTLE-W1-1-012418 (320-35355-1), FAY-D-384JERML-W1-1-012418 (320-35355-2), FAY-D-2528INDIA-W1-2-012418 (320-35355-3), FAY-D-661PINEA-W1-1-012418 (320-35355-4), FAY-D-701PINEA-W1-1-012418 (320-35355-5), FAY-D-701PINEA-W1-2-012418 (320-35355-6), FAY-D-663PINEA-W1-1-012418 (320-35355-7), FAY-D-7166BUTLE-W2-1-012418 (320-35355-8), FAY-D-7166BUTLE-W1-1-012418 (320-35355-9) and FAY-D-7166BUTLE-W1-1-012418D (320-35355-10) were analyzed for Perfluorinated Hydrocarbons in accordance with DV-LC-0012. The samples were prepared on 01/29/2018 and analyzed on 01/30/2018.

Calibration 9 (STD125) has been included in the raw data, but was not used in the Initial Calibration (ICAL).

Reporting limits have been adjusted accordingly for the initial volumes extracted.

No analytical or quality issues were noted, other than those described in the Definitions/Glossary page.

Fluorochemical Characterization of Water Samples Analytical Results

Chemours Sample Identification	TestAmerica Sample Identification	Collection Date/Time	Date Sample Received by TestAmerica	Analysis Date	HFPO-DA# (ug/L**)
FAY-D-7154BUTLE-W1-1-012418	320-35355-1	1/24/2018 15:19	1/25/2018	1/30/2018	0.047

HFPO-DA – hexafluoropropylene oxide dimer acid, analyzed by method DV-LC-0012, Revision 14.

< = less than the stated value

** ug/L – micrograms/liter (parts per billion)

DEFINITIONS:

Reporting Limit (RL) for the procedure is approximately 0.010ug/L.

RESULTS ARE CALCULATED ACCORDING TO THE FOLLOWING CRITERIA:

For samples analyzed in duplicate:

If the sample and laboratory duplicate are greater than 5X RL, the relative percent difference (RPD) is less than 20, the average value is reported. If the RPD is greater than 20, the higher value is reported.

If the sample or laboratory duplicate are less than 5X RL, and the absolute difference is less than RL, the average value is reported. If the absolute difference is greater than the RL, the higher value is reported.

Matrix Spike Recoveries:**Acceptable Range: 70%-130%**

TestAmerica Sample ID	Matrix Spike Recoveries
320-35355-9	91%

SUBMITTED BY:

2/8/2018

Carissa Cumine, Project Manager

Date

Fluorochemical Characterization of Water Samples Analytical Results

Chemours Sample Identification	TestAmerica Sample Identification	Collection Date/Time	Date Sample Received by TestAmerica	Analysis Date	HFPO-DA# (ug/L**)
FAY-D-384JERML-W1-1-012418	320-35355-2	1/24/2018 13:35	1/25/2018	1/30/2018	<0.010

HFPO-DA – hexafluoropropylene oxide dimer acid, analyzed by method DV-LC-0012, Revision 14.

< = less than the stated value

** ug/L – micrograms/liter (parts per billion)

DEFINITIONS:

Reporting Limit (RL) for the procedure is approximately 0.010ug/L.

RESULTS ARE CALCULATED ACCORDING TO THE FOLLOWING CRITERIA:

For samples analyzed in duplicate:

If the sample and laboratory duplicate are greater than 5X RL, the relative percent difference (RPD) is less than 20, the average value is reported. If the RPD is greater than 20, the higher value is reported.

If the sample or laboratory duplicate are less than 5X RL, and the absolute difference is less than RL, the average value is reported. If the absolute difference is greater than the RL, the higher value is reported.

Matrix Spike Recoveries:**Acceptable Range: 70%-130%**

TestAmerica Sample ID	Matrix Spike Recoveries
320-35355-9	91%

SUBMITTED BY:

2/8/2018

Carissa Cumine, Project Manager

Date

Fluorochemical Characterization of Water Samples Analytical Results

Chemours Sample Identification	TestAmerica Sample Identification	Collection Date/Time	Date Sample Received by TestAmerica	Analysis Date	HFPO-DA# (ug/L**)
FAY-D-2528INDIA-W1-2-012418	320-35355-3	1/24/2018 8:56	1/25/2018	1/30/2018	<0.010

HFPO-DA – hexafluoropropylene oxide dimer acid, analyzed by method DV-LC-0012, Revision 14.

< = less than the stated value

** ug/L – micrograms/liter (parts per billion)

DEFINITIONS:

Reporting Limit (RL) for the procedure is approximately 0.010ug/L.

RESULTS ARE CALCULATED ACCORDING TO THE FOLLOWING CRITERIA:

For samples analyzed in duplicate:

If the sample and laboratory duplicate are greater than 5X RL, the relative percent difference (RPD) is less than 20, the average value is reported. If the RPD is greater than 20, the higher value is reported.

If the sample or laboratory duplicate are less than 5X RL, and the absolute difference is less than RL, the average value is reported. If the absolute difference is greater than the RL, the higher value is reported.

Matrix Spike Recoveries:**Acceptable Range: 70%-130%**

TestAmerica Sample ID	Matrix Spike Recoveries
320-35355-9	91%

SUBMITTED BY:

2/8/2018

Carissa Cumine, Project Manager

Date

Fluorochemical Characterization of Water Samples Analytical Results

Chemours Sample Identification	TestAmerica Sample Identification	Collection Date/Time	Date Sample Received by TestAmerica	Analysis Date	HFPO-DA# (ug/L**)
FAY-D-661PINEA-W1-1-012418	320-35355-4	1/24/2018 9:40	1/25/2018	1/30/2018	0.040

HFPO-DA – hexafluoropropylene oxide dimer acid, analyzed by method DV-LC-0012, Revision 14.

< = less than the stated value

** ug/L – micrograms/liter (parts per billion)

DEFINITIONS:

Reporting Limit (RL) for the procedure is approximately 0.010ug/L.

RESULTS ARE CALCULATED ACCORDING TO THE FOLLOWING CRITERIA:


For samples analyzed in duplicate:

If the sample and laboratory duplicate are greater than 5X RL, the relative percent difference (RPD) is less than 20, the average value is reported. If the RPD is greater than 20, the higher value is reported.

If the sample or laboratory duplicate are less than 5X RL, and the absolute difference is less than RL, the average value is reported. If the absolute difference is greater than the RL, the higher value is reported.

Matrix Spike Recoveries:**Acceptable Range: 70%-130%**

TestAmerica Sample ID	Matrix Spike Recoveries
320-35355-9	91%

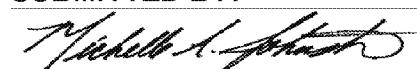SUBMITTED BY:

3/29/2018

Michelle Johnston, Project Manager

Date

Fluorochemical Characterization of Water Samples Analytical Results

Chemours Sample Identification	TestAmerica Sample Identification	Collection Date/Time	Date Sample Received by TestAmerica	Analysis Date	HFPO-DA# (ug/L**)
FAY-D-701PINEA-W1-1-012418	320-35355-5	1/24/2018 9:58	1/25/2018	1/30/2018	<0.010

HFPO-DA – hexafluoropropylene oxide dimer acid, analyzed by method DV-LC-0012, Revision 14.

< = less than the stated value

** ug/L – micrograms/liter (parts per billion)

DEFINITIONS:

Reporting Limit (RL) for the procedure is approximately 0.010ug/L.

RESULTS ARE CALCULATED ACCORDING TO THE FOLLOWING CRITERIA:

For samples analyzed in duplicate:

If the sample and laboratory duplicate are greater than 5X RL, the relative percent difference (RPD) is less than 20, the average value is reported. If the RPD is greater than 20, the higher value is reported.

If the sample or laboratory duplicate are less than 5X RL, and the absolute difference is less than RL, the average value is reported. If the absolute difference is greater than the RL, the higher value is reported.

Matrix Spike Recoveries:

Acceptable Range: 70%-130%

TestAmerica Sample ID	Matrix Spike Recoveries
320-35355-9	91%

SUBMITTED BY:

2/8/2018

Carissa Cumine, Project Manager

Date

Fluorochemical Characterization of Water Samples Analytical Results

Chemours Sample Identification	TestAmerica Sample Identification	Collection Date/Time	Date Sample Received by TestAmerica	Analysis Date	HFPO-DA# (ug/L**)
FAY-D-701PINEA-W1-2-012418	320-35355-6	1/24/2018 10:05	1/25/2018	1/30/2018	<0.010

HFPO-DA – hexafluoropropylene oxide dimer acid, analyzed by method DV-LC-0012, Revision 14.

< = less than the stated value

** ug/L – micrograms/liter (parts per billion)

DEFINITIONS:

Reporting Limit (RL) for the procedure is approximately 0.010ug/L.

RESULTS ARE CALCULATED ACCORDING TO THE FOLLOWING CRITERIA:

For samples analyzed in duplicate:

If the sample and laboratory duplicate are greater than 5X RL, the relative percent difference (RPD) is less than 20, the average value is reported. If the RPD is greater than 20, the higher value is reported.

If the sample or laboratory duplicate are less than 5X RL, and the absolute difference is less than RL, the average value is reported. If the absolute difference is greater than the RL, the higher value is reported.

Matrix Spike Recoveries:**Acceptable Range: 70%-130%**

TestAmerica Sample ID	Matrix Spike Recoveries
320-35355-9	91%

SUBMITTED BY:

2/8/2018

Carissa Cumine, Project Manager

Date

Fluorochemical Characterization of Water Samples Analytical Results

Chemours Sample Identification	TestAmerica Sample Identification	Collection Date/Time	Date Sample Received by TestAmerica	Analysis Date	HFPO-DA# (ug/L**)
FAY-D-663PINEA-W1-1-012418	320-35355-7	1/24/2018 10:29	1/25/2018	1/30/2018	0.015

HFPO-DA – hexafluoropropylene oxide dimer acid, analyzed by method DV-LC-0012, Revision 14.

< = less than the stated value

** ug/L – micrograms/liter (parts per billion)

DEFINITIONS:

Reporting Limit (RL) for the procedure is approximately 0.010ug/L.

RESULTS ARE CALCULATED ACCORDING TO THE FOLLOWING CRITERIA:

For samples analyzed in duplicate:

If the sample and laboratory duplicate are greater than 5X RL, the relative percent difference (RPD) is less than 20, the average value is reported. If the RPD is greater than 20, the higher value is reported.

If the sample or laboratory duplicate are less than 5X RL, and the absolute difference is less than RL, the average value is reported. If the absolute difference is greater than the RL, the higher value is reported.

Matrix Spike Recoveries:

Acceptable Range: 70%-130%

TestAmerica Sample ID	Matrix Spike Recoveries
320-35355-9	91%

SUBMITTED BY:

2/8/2018

Carissa Cumine, Project Manager

Date

Fluorochemical Characterization of Water Samples Analytical Results

Chemours Sample Identification	TestAmerica Sample Identification	Collection Date/Time	Date Sample Received by TestAmerica	Analysis Date	HFPO-DA# (ug/L**)
FAY-D-7166BUTLE-W2-1-012418	320-35355-8	1/24/2018 15:43	1/25/2018	1/30/2018	0.021

HFPO-DA – hexafluoropropylene oxide dimer acid, analyzed by method DV-LC-0012, Revision 14.

< = less than the stated value

** ug/L – micrograms/liter (parts per billion)

DEFINITIONS:

Reporting Limit (RL) for the procedure is approximately 0.010ug/L.

RESULTS ARE CALCULATED ACCORDING TO THE FOLLOWING CRITERIA:

For samples analyzed in duplicate:

If the sample and laboratory duplicate are greater than 5X RL, the relative percent difference (RPD) is less than 20, the average value is reported. If the RPD is greater than 20, the higher value is reported.

If the sample or laboratory duplicate are less than 5X RL, and the absolute difference is less than RL, the average value is reported. If the absolute difference is greater than the RL, the higher value is reported.

Matrix Spike Recoveries:**Acceptable Range: 70%-130%**

TestAmerica Sample ID	Matrix Spike Recoveries
320-35355-9	91%

SUBMITTED BY:

2/8/2018

Carissa Cumine, Project Manager

Date

Fluorochemical Characterization of Water Samples Analytical Results

Chemours Sample Identification	TestAmerica Sample Identification	Collection Date/Time	Date Sample Received by TestAmerica	Analysis Date	HFPO-DA# (ug/L**)
FAY-D-7166BUTLE-W1-1-012418	320-35355-9	1/24/2018 15:35	1/25/2018	1/30/2018	0.017

HFPO-DA – hexafluoropropylene oxide dimer acid, analyzed by method DV-LC-0012, Revision 14.

< = less than the stated value

** ug/L – micrograms/liter (parts per billion)

DEFINITIONS:

Reporting Limit (RL) for the procedure is approximately 0.010ug/L.

RESULTS ARE CALCULATED ACCORDING TO THE FOLLOWING CRITERIA:

For samples analyzed in duplicate:

If the sample and laboratory duplicate are greater than 5X RL, the relative percent difference (RPD) is less than 20, the average value is reported. If the RPD is greater than 20, the higher value is reported.

If the sample or laboratory duplicate are less than 5X RL, and the absolute difference is less than RL, the average value is reported. If the absolute difference is greater than the RL, the higher value is reported.

Matrix Spike Recoveries:**Acceptable Range: 70%-130%**

TestAmerica Sample ID	Matrix Spike Recoveries
320-35355-9	91%

SUBMITTED BY:

Cumine
Carissa Cumine, Project Manager

2/8/2018

Date

Fluorochemical Characterization of Water Samples Analytical Results

Chemours Sample Identification	TestAmerica Sample Identification	Collection Date/Time	Date Sample Received by TestAmerica	Analysis Date	HFPO-DA# (ug/L**)
FAY-D-7166BUTLE-W1-1-012418D	320-35355-10	1/24/2018 15:35	1/25/2018	1/30/2018	0.021

HFPO-DA – hexafluoropropylene oxide dimer acid, analyzed by method DV-LC-0012, Revision 14.

< = less than the stated value

** ug/L – micrograms/liter (parts per billion)

DEFINITIONS:

Reporting Limit (RL) for the procedure is approximately 0.010ug/L.

RESULTS ARE CALCULATED ACCORDING TO THE FOLLOWING CRITERIA:

For samples analyzed in duplicate:

If the sample and laboratory duplicate are greater than 5X RL, the relative percent difference (RPD) is less than 20, the average value is reported. If the RPD is greater than 20, the higher value is reported.

If the sample or laboratory duplicate are less than 5X RL, and the absolute difference is less than RL, the average value is reported. If the absolute difference is greater than the RL, the higher value is reported.

Matrix Spike Recoveries:**Acceptable Range: 70%-130%**

TestAmerica Sample ID	Matrix Spike Recoveries
320-35355-9	91%

SUBMITTED BY:

2/8/2018

Carissa Cumine, Project Manager

Date

Executive Summary

Client: Chemours Company FC, LLC The

Job Number: 320-35355-1

8321A : HFPO-DA

Lab Sample ID	Client Sample ID	Analyte	Individual Result (ug/L)	Final Result (ug/L)	RL
320-35355-1	FAY-D-7154BUTLE-W1-1-012418	HFPO-DA	0.047	0.047	0.010
320-35355-2	FAY-D-384JERML-W1-1-012418	HFPO-DA	<0.010	<0.010	0.010
320-35355-3	FAY-D-2528INDIA-W1-2-012418	HFPO-DA	<0.010	<0.010	0.010
320-35355-4	FAY-D-661PINEA-W1-1-012418	HFPO-DA	0.040	0.040	0.010
320-35355-5	FAY-D-701PINEA-W1-1-012418	HFPO-DA	<0.010	<0.010	0.010
320-35355-6	FAY-D-701PINEA-W1-2-012418	HFPO-DA	<0.010	<0.010	0.010
320-35355-7	FAY-D-663PINEA-W1-1-012418	HFPO-DA	0.015	0.015	0.010
320-35355-8	FAY-D-7166BUTLE-W2-1-012418	HFPO-DA	0.021	0.021	0.010
320-35355-9	FAY-D-7166BUTLE-W1-1-012418	HFPO-DA	0.017	0.017	0.010
320-35355-9 DU	FAY-D-7166BUTLE-W1-1-012418	HFPO-DA	0.017		0.010
320-35355-10	FAY-D-7166BUTLE-W1-1-012418D	HFPO-DA	0.021	0.021	0.010

(a) Method 8321A

(b) DUP or REP indicates a laboratory duplicate.

(c) If the sample and laboratory duplicate are both greater than 5X the RL and the relative percent difference (RPD) is less than 20, the average value is reported. If the RPD is greater than 20, the higher of the sample and laboratory duplicate value is reported. If the sample and/or laboratory duplicate are less than 5X the RL, and the absolute difference between the sample and laboratory duplicate is less than the RL, the average value is reported. If the absolute difference is greater than the RL, the higher of the sample and laboratory duplicate value is reported. If either the sample or the duplicate result is greater than or equal to the RL and the other is less than the RL, then the higher of the two is reported.

(d) Moisture Determined by ASTM D2216.

(e) Reporting Limit (RL) = The concentration equivalent to the low calibration standard.

Detection Summary

Client: Chemours Company FC, LLC The
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 320-35355-1

Client Sample ID: FAY-D-7154BUTLE-W1-1-012418

Lab Sample ID: 320-35355-1

Analyte	Result	Qualifier	RL	MDL	Unit	Dil Fac	D	Method	Prep Type
HFPO-DA	0.047		0.010		ug/L	1		8321A	Total/NA

Client Sample ID: FAY-D-384JERML-W1-1-012418

Lab Sample ID: 320-35355-2

No Detections.

Client Sample ID: FAY-D-2528INDIA-W1-2-012418

Lab Sample ID: 320-35355-3

No Detections.

Client Sample ID: FAY-D-661PINEA-W1-1-012418

Lab Sample ID: 320-35355-4

Analyte	Result	Qualifier	RL	MDL	Unit	Dil Fac	D	Method	Prep Type
HFPO-DA	0.040		0.010		ug/L	1		8321A	Total/NA

Client Sample ID: FAY-D-701PINEA-W1-1-012418

Lab Sample ID: 320-35355-5

No Detections.

Client Sample ID: FAY-D-701PINEA-W1-2-012418

Lab Sample ID: 320-35355-6

No Detections.

Client Sample ID: FAY-D-663PINEA-W1-1-012418

Lab Sample ID: 320-35355-7

Analyte	Result	Qualifier	RL	MDL	Unit	Dil Fac	D	Method	Prep Type
HFPO-DA	0.015		0.010		ug/L	1		8321A	Total/NA

Client Sample ID: FAY-D-7166BUTLE-W2-1-012418

Lab Sample ID: 320-35355-8

Analyte	Result	Qualifier	RL	MDL	Unit	Dil Fac	D	Method	Prep Type
HFPO-DA	0.021		0.010		ug/L	1		8321A	Total/NA

Client Sample ID: FAY-D-7166BUTLE-W1-1-012418

Lab Sample ID: 320-35355-9

Analyte	Result	Qualifier	RL	MDL	Unit	Dil Fac	D	Method	Prep Type
HFPO-DA	0.017		0.010		ug/L	1		8321A	Total/NA

Client Sample ID: FAY-D-7166BUTLE-W1-1-012418D

Lab Sample ID: 320-35355-10

Analyte	Result	Qualifier	RL	MDL	Unit	Dil Fac	D	Method	Prep Type
HFPO-DA	0.021		0.010		ug/L	1		8321A	Total/NA

This Detection Summary does not include radiochemical test results.

TestAmerica Sacramento

Client Sample Results

Client: Chemours Company FC, LLC The
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 320-35355-1

Client Sample ID: FAY-D-7154BUTLE-W1-1-012418

Lab Sample ID: 320-35355-1

Date Collected: 01/24/18 15:19

Matrix: Water

Date Received: 01/25/18 09:30

Method: 8321A - HFPO-DA

Analyte	Result	Qualifier	RL	MDL	Unit	D	Prepared	Analyzed	Dil Fac
HFPO-DA	0.047		0.010		ug/L		01/29/18 21:37	01/30/18 09:28	1
Surrogate	%Recovery	Qualifier	Limits				Prepared	Analyzed	Dil Fac
13C3 HFPO-DA	103		50 - 200				01/29/18 21:37	01/30/18 09:28	1

Client Sample ID: FAY-D-384JERML-W1-1-012418

Lab Sample ID: 320-35355-2

Date Collected: 01/24/18 13:35

Matrix: Water

Date Received: 01/25/18 09:30

Method: 8321A - HFPO-DA

Analyte	Result	Qualifier	RL	MDL	Unit	D	Prepared	Analyzed	Dil Fac
HFPO-DA	<0.010		0.010		ug/L		01/29/18 21:37	01/30/18 09:32	1
Surrogate	%Recovery	Qualifier	Limits				Prepared	Analyzed	Dil Fac
13C3 HFPO-DA	107		50 - 200				01/29/18 21:37	01/30/18 09:32	1

Client Sample ID: FAY-D-2528INDIA-W1-2-012418

Lab Sample ID: 320-35355-3

Date Collected: 01/24/18 08:56

Matrix: Water

Date Received: 01/25/18 09:30

Method: 8321A - HFPO-DA

Analyte	Result	Qualifier	RL	MDL	Unit	D	Prepared	Analyzed	Dil Fac
HFPO-DA	<0.010		0.010		ug/L		01/29/18 21:37	01/30/18 09:35	1
Surrogate	%Recovery	Qualifier	Limits				Prepared	Analyzed	Dil Fac
13C3 HFPO-DA	104		50 - 200				01/29/18 21:37	01/30/18 09:35	1

Client Sample ID: FAY-D-661PINEA-W1-1-012418

Lab Sample ID: 320-35355-4

Date Collected: 01/24/18 09:40

Matrix: Water

Date Received: 01/25/18 09:30

Method: 8321A - HFPO-DA

Analyte	Result	Qualifier	RL	MDL	Unit	D	Prepared	Analyzed	Dil Fac
HFPO-DA	0.040		0.010		ug/L		01/29/18 21:37	01/30/18 09:38	1
Surrogate	%Recovery	Qualifier	Limits				Prepared	Analyzed	Dil Fac
13C3 HFPO-DA	97		50 - 200				01/29/18 21:37	01/30/18 09:38	1

Client Sample ID: FAY-D-701PINEA-W1-1-012418

Lab Sample ID: 320-35355-5

Date Collected: 01/24/18 09:58

Matrix: Water

Date Received: 01/25/18 09:30

Method: 8321A - HFPO-DA

Analyte	Result	Qualifier	RL	MDL	Unit	D	Prepared	Analyzed	Dil Fac
HFPO-DA	<0.010		0.010		ug/L		01/29/18 21:37	01/30/18 09:45	1
Surrogate	%Recovery	Qualifier	Limits				Prepared	Analyzed	Dil Fac
13C3 HFPO-DA	102		50 - 200				01/29/18 21:37	01/30/18 09:45	1

TestAmerica Sacramento

Client Sample Results

Client: Chemours Company FC, LLC The
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 320-35355-1

Client Sample ID: FAY-D-701PINEA-W1-2-012418

Lab Sample ID: 320-35355-6

Matrix: Water

Date Collected: 01/24/18 10:05

Date Received: 01/25/18 09:30

Method: 8321A - HFPO-DA

Analyte	Result	Qualifier	RL	MDL	Unit	D	Prepared	Analyzed	Dil Fac
HFPO-DA	<0.010		0.010		ug/L	D	01/29/18 21:37	01/30/18 09:48	1
<i>Surrogate</i>	<i>%Recovery</i>	<i>Qualifier</i>	<i>Limits</i>				<i>Prepared</i>	<i>Analyzed</i>	<i>Dil Fac</i>
13C3 HFPO-DA	125		50 - 200				01/29/18 21:37	01/30/18 09:48	1

Client Sample ID: FAY-D-663PINEA-W1-1-012418

Lab Sample ID: 320-35355-7

Matrix: Water

Date Collected: 01/24/18 10:29

Date Received: 01/25/18 09:30

Method: 8321A - HFPO-DA

Analyte	Result	Qualifier	RL	MDL	Unit	D	Prepared	Analyzed	Dil Fac
HFPO-DA	0.015		0.010		ug/L	D	01/29/18 21:37	01/30/18 09:51	1
<i>Surrogate</i>	<i>%Recovery</i>	<i>Qualifier</i>	<i>Limits</i>				<i>Prepared</i>	<i>Analyzed</i>	<i>Dil Fac</i>
13C3 HFPO-DA	123		50 - 200				01/29/18 21:37	01/30/18 09:51	1

Client Sample ID: FAY-D-7166BUTLE-W2-1-012418

Lab Sample ID: 320-35355-8

Matrix: Water

Date Collected: 01/24/18 15:43

Date Received: 01/25/18 09:30

Method: 8321A - HFPO-DA

Analyte	Result	Qualifier	RL	MDL	Unit	D	Prepared	Analyzed	Dil Fac
HFPO-DA	0.021		0.010		ug/L	D	01/29/18 21:37	01/30/18 09:54	1
<i>Surrogate</i>	<i>%Recovery</i>	<i>Qualifier</i>	<i>Limits</i>				<i>Prepared</i>	<i>Analyzed</i>	<i>Dil Fac</i>
13C3 HFPO-DA	100		50 - 200				01/29/18 21:37	01/30/18 09:54	1

Client Sample ID: FAY-D-7166BUTLE-W1-1-012418

Lab Sample ID: 320-35355-9

Matrix: Water

Date Collected: 01/24/18 15:35

Date Received: 01/25/18 09:30

Method: 8321A - HFPO-DA

Analyte	Result	Qualifier	RL	MDL	Unit	D	Prepared	Analyzed	Dil Fac
HFPO-DA	0.017		0.010		ug/L	D	01/29/18 21:37	01/30/18 09:58	1
<i>Surrogate</i>	<i>%Recovery</i>	<i>Qualifier</i>	<i>Limits</i>				<i>Prepared</i>	<i>Analyzed</i>	<i>Dil Fac</i>
13C3 HFPO-DA	90		50 - 200				01/29/18 21:37	01/30/18 09:58	1

Client Sample ID: FAY-D-7166BUTLE-W1-1-012418D

Lab Sample ID: 320-35355-10

Matrix: Water

Date Collected: 01/24/18 15:35

Date Received: 01/25/18 09:30

Method: 8321A - HFPO-DA

Analyte	Result	Qualifier	RL	MDL	Unit	D	Prepared	Analyzed	Dil Fac
HFPO-DA	0.021		0.010		ug/L	D	01/29/18 21:37	01/30/18 10:07	1
<i>Surrogate</i>	<i>%Recovery</i>	<i>Qualifier</i>	<i>Limits</i>				<i>Prepared</i>	<i>Analyzed</i>	<i>Dil Fac</i>
13C3 HFPO-DA	103		50 - 200				01/29/18 21:37	01/30/18 10:07	1

TestAmerica Sacramento

Default Detection Limits

Client: Chemours Company FC, LLC The
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 320-35355-1

Method: 8321A - HFPO-DA

Prep: 3535

Analyte	RL	MDL	Units	Method
HFPO-DA	0.010	0.0051	ug/L	8321A

Surrogate Summary

Client: Chemours Company FC, LLC The
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 320-35355-1

Method: 8321A - HFPO-DA

Matrix: Water

Prep Type: Total/NA

Percent Surrogate Recovery (Acceptance Limits)

HFPODA

(50-200)

Lab Sample ID	Client Sample ID	HFPODA (50-200)
320-35355-1	FAY-D-7154BUTLE-W1-1-0124	103
320-35355-2	FAY-D-384JERML-W1-1-01241	107
320-35355-3	FAY-D-2528INDIA-W1-2-01241	104
320-35355-4	FAY-D-661PINEA-W1-1-012418	97
320-35355-5	FAY-D-701PINEA-W1-1-012418	102
320-35355-6	FAY-D-701PINEA-W1-2-012418	125
320-35355-7	FAY-D-663PINEA-W1-1-012418	123
320-35355-8	FAY-D-7166BUTLE-W2-1-0124	100
320-35355-9	FAY-D-7166BUTLE-W1-1-0124	90
320-35355-9 DU	FAY-D-7166BUTLE-W1-1-0124	100
320-35355-9 MS	FAY-D-7166BUTLE-W1-1-0124	123
320-35355-10	FAY-D-7166BUTLE-W1-1-0124	103
DLCK 280-390728/12	Lab Control Sample	102
LCS 280-403139/2-A	Lab Control Sample	107
LCSD 280-403139/3-A	Lab Control Sample Dup	105
LLCS 280-403139/4-A	Lab Control Sample	128
MB 280-403139/1-A	Method Blank	111

Surrogate Legend

HFPODA = 13C3 HFPO-DA

TestAmerica Sacramento

QC Sample Results

Client: Chemours Company FC, LLC The
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 320-35355-1

Method: 8321A - HFPO-DA

Lab Sample ID: DLCK 280-390728/12

Matrix: Water

Analysis Batch: 390728

Client Sample ID: Lab Control Sample
Prep Type: Total/NA

Analyte	Spike Added	DLCK Result	DLCK Qualifier	Unit	D	% Rec.	% Rec. Limits
HFPO-DA	0.250	<0.50		ug/L		78	70 - 130
Surrogate	DLCK %Recovery	DLCK Qualifier	Limits				
13C3 HFPO-DA	102		50 - 200				

Lab Sample ID: MB 280-403139/1-A

Matrix: Water

Analysis Batch: 403246

Client Sample ID: Method Blank
Prep Type: Total/NA
Prep Batch: 403139

Analyte	MB Result	MB Qualifier	RL	MDL	Unit	D	Prepared	Analyzed	Dil Fac
HFPO-DA	<0.010		0.010		ug/L		01/29/18 21:37	01/30/18 09:15	1
Surrogate	MB %Recovery	MB Qualifier	Limits				Prepared	Analyzed	Dil Fac
13C3 HFPO-DA	111		50 - 200				01/29/18 21:37	01/30/18 09:15	1

Lab Sample ID: LCS 280-403139/2-A

Matrix: Water

Analysis Batch: 403246

Client Sample ID: Lab Control Sample
Prep Type: Total/NA
Prep Batch: 403139

Analyte	Spike Added	LCS Result	LCS Qualifier	Unit	D	% Rec.	% Rec. Limits
HFPO-DA	0.200	0.225		ug/L		113	70 - 130
Surrogate	LCS %Recovery	LCS Qualifier	Limits				
13C3 HFPO-DA	107		50 - 200				

Lab Sample ID: LCSD 280-403139/3-A

Matrix: Water

Analysis Batch: 403246

Client Sample ID: Lab Control Sample Dup
Prep Type: Total/NA
Prep Batch: 403139

Analyte	Spike Added	LCSD Result	LCSD Qualifier	Unit	D	% Rec.	% Rec. Limits	RPD	Limit
HFPO-DA	0.200	0.207		ug/L		103	70 - 130	8	20
Surrogate	LCSD %Recovery	LCSD Qualifier	Limits						
13C3 HFPO-DA	105		50 - 200						

Lab Sample ID: LLCS 280-403139/4-A

Matrix: Water

Analysis Batch: 403246

Client Sample ID: Lab Control Sample
Prep Type: Total/NA
Prep Batch: 403139

Analyte	Spike Added	LLCS Result	LLCS Qualifier	Unit	D	% Rec.	% Rec. Limits
HFPO-DA	0.0200	0.0147		ug/L		74	70 - 130
Surrogate	LLCS %Recovery	LLCS Qualifier	Limits				
13C3 HFPO-DA	128		50 - 200				

TestAmerica Sacramento

QC Sample Results

Client: Chemours Company FC, LLC The
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 320-35355-1

Method: 8321A - HFPO-DA (Continued)

Lab Sample ID: 320-35355-9 MS

Matrix: Water

Analysis Batch: 403246

Client Sample ID: FAY-D-7166BUTLE-W1-1-012418

Prep Type: Total/NA

Prep Batch: 403139

% Rec.

Limits

Analyte	Sample Result	Sample Qualifier	Spike Added	MS Result	MS Qualifier	Unit	D	% Rec.	RPD
HFPO-DA	0.017		0.203	0.202		ug/L	91	70 - 130	
Surrogate	MS %Recovery	MS Qualifier		MS Result	MS Qualifier	Unit	D	% Rec.	RPD
13C3 HFPO-DA	123			50 - 200					

Lab Sample ID: 320-35355-9 DU

Matrix: Water

Analysis Batch: 403246

Client Sample ID: FAY-D-7166BUTLE-W1-1-012418

Prep Type: Total/NA

Prep Batch: 403139

RPD

Limit

Analyte	Sample Result	Sample Qualifier		DU Result	DU Qualifier	Unit	D	RPD	Limit
HFPO-DA	0.017			0.0171		ug/L		1	20
Surrogate	DU %Recovery	DU Qualifier		DU Result	DU Qualifier	Unit	D	RPD	Limit
13C3 HFPO-DA	100			50 - 200					

QC Association Summary

Client: Chemours Company FC, LLC The
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 320-35355-1

LCMS

Analysis Batch: 390728

Lab Sample ID	Client Sample ID	Prep Type	Matrix	Method	Prep Batch
DLCK 280-390728/12	Lab Control Sample	Total/NA	Water	8321A	

Prep Batch: 403139

Lab Sample ID	Client Sample ID	Prep Type	Matrix	Method	Prep Batch
320-35355-1	FAY-D-7154BUTLE-W1-1-012418	Total/NA	Water	3535	
320-35355-2	FAY-D-384JERML-W1-1-012418	Total/NA	Water	3535	
320-35355-3	FAY-D-2528INDIA-W1-2-012418	Total/NA	Water	3535	
320-35355-4	FAY-D-661PINEA-W1-1-012418	Total/NA	Water	3535	
320-35355-5	FAY-D-701PINEA-W1-1-012418	Total/NA	Water	3535	
320-35355-6	FAY-D-701PINEA-W1-2-012418	Total/NA	Water	3535	
320-35355-7	FAY-D-663PINEA-W1-1-012418	Total/NA	Water	3535	
320-35355-8	FAY-D-7166BUTLE-W2-1-012418	Total/NA	Water	3535	
320-35355-9	FAY-D-7166BUTLE-W1-1-012418	Total/NA	Water	3535	
320-35355-10	FAY-D-7166BUTLE-W1-1-012418D	Total/NA	Water	3535	
MB 280-403139/1-A	Method Blank	Total/NA	Water	3535	
LCS 280-403139/2-A	Lab Control Sample	Total/NA	Water	3535	
LCSD 280-403139/3-A	Lab Control Sample Dup	Total/NA	Water	3535	
LLCS 280-403139/4-A	Lab Control Sample	Total/NA	Water	3535	
320-35355-9 MS	FAY-D-7166BUTLE-W1-1-012418	Total/NA	Water	3535	
320-35355-9 DU	FAY-D-7166BUTLE-W1-1-012418	Total/NA	Water	3535	

Analysis Batch: 403246

Lab Sample ID	Client Sample ID	Prep Type	Matrix	Method	Prep Batch
320-35355-1	FAY-D-7154BUTLE-W1-1-012418	Total/NA	Water	8321A	403139
320-35355-2	FAY-D-384JERML-W1-1-012418	Total/NA	Water	8321A	403139
320-35355-3	FAY-D-2528INDIA-W1-2-012418	Total/NA	Water	8321A	403139
320-35355-4	FAY-D-661PINEA-W1-1-012418	Total/NA	Water	8321A	403139
320-35355-5	FAY-D-701PINEA-W1-1-012418	Total/NA	Water	8321A	403139
320-35355-6	FAY-D-701PINEA-W1-2-012418	Total/NA	Water	8321A	403139
320-35355-7	FAY-D-663PINEA-W1-1-012418	Total/NA	Water	8321A	403139
320-35355-8	FAY-D-7166BUTLE-W2-1-012418	Total/NA	Water	8321A	403139
320-35355-9	FAY-D-7166BUTLE-W1-1-012418	Total/NA	Water	8321A	403139
320-35355-10	FAY-D-7166BUTLE-W1-1-012418D	Total/NA	Water	8321A	403139
MB 280-403139/1-A	Method Blank	Total/NA	Water	8321A	403139
LCS 280-403139/2-A	Lab Control Sample	Total/NA	Water	8321A	403139
LCSD 280-403139/3-A	Lab Control Sample Dup	Total/NA	Water	8321A	403139
LLCS 280-403139/4-A	Lab Control Sample	Total/NA	Water	8321A	403139
320-35355-9 MS	FAY-D-7166BUTLE-W1-1-012418	Total/NA	Water	8321A	403139
320-35355-9 DU	FAY-D-7166BUTLE-W1-1-012418	Total/NA	Water	8321A	403139

TestAmerica Sacramento

Lab Chronicle

Client: Chemours Company FC, LLC The
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 320-35355-1

Client Sample ID: FAY-D-7154BUTLE-W1-1-012418

Lab Sample ID: 320-35355-1

Date Collected: 01/24/18 15:19

Matrix: Water

Date Received: 01/25/18 09:30

Prep Type	Batch Type	Batch Method	Run	Dilution Factor	Batch Number	Prepared or Analyzed	Analyst	Lab
Total/NA	Prep	3535			403139	01/29/18 21:37	CDC	TAL DEN
Total/NA	Analysis	8321A		1	403246	01/30/18 09:28	AGCM	TAL DEN

Client Sample ID: FAY-D-384JERML-W1-1-012418

Lab Sample ID: 320-35355-2

Date Collected: 01/24/18 13:35

Matrix: Water

Date Received: 01/25/18 09:30

Prep Type	Batch Type	Batch Method	Run	Dilution Factor	Batch Number	Prepared or Analyzed	Analyst	Lab
Total/NA	Prep	3535			403139	01/29/18 21:37	CDC	TAL DEN
Total/NA	Analysis	8321A		1	403246	01/30/18 09:32	AGCM	TAL DEN

Client Sample ID: FAY-D-2528INDIA-W1-2-012418

Lab Sample ID: 320-35355-3

Date Collected: 01/24/18 08:56

Matrix: Water

Date Received: 01/25/18 09:30

Prep Type	Batch Type	Batch Method	Run	Dilution Factor	Batch Number	Prepared or Analyzed	Analyst	Lab
Total/NA	Prep	3535			403139	01/29/18 21:37	CDC	TAL DEN
Total/NA	Analysis	8321A		1	403246	01/30/18 09:35	AGCM	TAL DEN

Client Sample ID: FAY-D-661PINEA-W1-1-012418

Lab Sample ID: 320-35355-4

Date Collected: 01/24/18 09:40

Matrix: Water

Date Received: 01/25/18 09:30

Prep Type	Batch Type	Batch Method	Run	Dilution Factor	Batch Number	Prepared or Analyzed	Analyst	Lab
Total/NA	Prep	3535			403139	01/29/18 21:37	CDC	TAL DEN
Total/NA	Analysis	8321A		1	403246	01/30/18 09:38	AGCM	TAL DEN

Client Sample ID: FAY-D-701PINEA-W1-1-012418

Lab Sample ID: 320-35355-5

Date Collected: 01/24/18 09:58

Matrix: Water

Date Received: 01/25/18 09:30

Prep Type	Batch Type	Batch Method	Run	Dilution Factor	Batch Number	Prepared or Analyzed	Analyst	Lab
Total/NA	Prep	3535			403139	01/29/18 21:37	CDC	TAL DEN
Total/NA	Analysis	8321A		1	403246	01/30/18 09:45	AGCM	TAL DEN

Client Sample ID: FAY-D-701PINEA-W1-2-012418

Lab Sample ID: 320-35355-6

Date Collected: 01/24/18 10:05

Matrix: Water

Date Received: 01/25/18 09:30

Prep Type	Batch Type	Batch Method	Run	Dilution Factor	Batch Number	Prepared or Analyzed	Analyst	Lab
Total/NA	Prep	3535			403139	01/29/18 21:37	CDC	TAL DEN
Total/NA	Analysis	8321A		1	403246	01/30/18 09:48	AGCM	TAL DEN

TestAmerica Sacramento

Lab Chronicle

Client: Chemours Company FC, LLC The
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 320-35355-1

Client Sample ID: FAY-D-663PINEA-W1-1-012418

Lab Sample ID: 320-35355-7

Date Collected: 01/24/18 10:29

Matrix: Water

Date Received: 01/25/18 09:30

Prep Type	Batch Type	Batch Method	Run	Dilution Factor	Batch Number	Prepared or Analyzed	Analyst	Lab
Total/NA	Prep	3535			403139	01/29/18 21:37	CDC	TAL DEN
Total/NA	Analysis	8321A		1	403246	01/30/18 09:51	AGCM	TAL DEN

Client Sample ID: FAY-D-7166BUTLE-W2-1-012418

Lab Sample ID: 320-35355-8

Date Collected: 01/24/18 15:43

Matrix: Water

Date Received: 01/25/18 09:30

Prep Type	Batch Type	Batch Method	Run	Dilution Factor	Batch Number	Prepared or Analyzed	Analyst	Lab
Total/NA	Prep	3535			403139	01/29/18 21:37	CDC	TAL DEN
Total/NA	Analysis	8321A		1	403246	01/30/18 09:54	AGCM	TAL DEN

Client Sample ID: FAY-D-7166BUTLE-W1-1-012418

Lab Sample ID: 320-35355-9

Date Collected: 01/24/18 15:35

Matrix: Water

Date Received: 01/25/18 09:30

Prep Type	Batch Type	Batch Method	Run	Dilution Factor	Batch Number	Prepared or Analyzed	Analyst	Lab
Total/NA	Prep	3535			403139	01/29/18 21:37	CDC	TAL DEN
Total/NA	Analysis	8321A		1	403246	01/30/18 09:58	AGCM	TAL DEN

Client Sample ID: FAY-D-7166BUTLE-W1-1-012418D

Lab Sample ID: 320-35355-10

Date Collected: 01/24/18 15:35

Matrix: Water

Date Received: 01/25/18 09:30

Prep Type	Batch Type	Batch Method	Run	Dilution Factor	Batch Number	Prepared or Analyzed	Analyst	Lab
Total/NA	Prep	3535			403139	01/29/18 21:37	CDC	TAL DEN
Total/NA	Analysis	8321A		1	403246	01/30/18 10:07	AGCM	TAL DEN

Laboratory References:

TAL DEN = TestAmerica Denver, 4955 Yarrow Street, Arvada, CO 80002, TEL (303)736-0100

TestAmerica Sacramento

Accreditation/Certification Summary

Client: Chemours Company FC, LLC The
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 320-35355-1

Laboratory: TestAmerica Sacramento

All accreditations/certifications held by this laboratory are listed. Not all accreditations/certifications are applicable to this report.

Authority	Program	EPA Region	Identification Number	Expiration Date
Alaska (UST)	State Program	10	17-020	01-20-21
Arizona	State Program	9	AZ0708	08-11-18
Arkansas DEQ	State Program	6	88-0691	06-17-18
California	State Program	9	2897	01-31-19
Colorado	State Program	8	CA00044	08-31-18
Connecticut	State Program	1	PH-0691	06-30-19
Florida	NELAP	4	E87570	06-30-18
Georgia	State Program	4	N/A	01-28-19
Hawaii	State Program	9	N/A	01-29-19
Illinois	NELAP	5	200060	03-17-18
Kansas	NELAP	7	E-10375	10-31-18
L-A-B	DoD ELAP		L2468	01-20-21
Louisiana	NELAP	6	30612	06-30-18
Maine	State Program	1	CA0004	04-14-18
Michigan	State Program	5	9947	01-31-18 *
Nevada	State Program	9	CA00044	07-31-18
New Hampshire	NELAP	1	2997	04-18-18
New Jersey	NELAP	2	CA005	06-30-18
New York	NELAP	2	11666	04-01-18
Oregon	NELAP	10	4040	01-29-20
Pennsylvania	NELAP	3	68-01272	03-31-18
Texas	NELAP	6	T104704399	05-31-18
US Fish & Wildlife	Federal		LE148388-0	07-31-18
USDA	Federal		P330-11-00436	01-17-21
USEPA UCMR	Federal	1	CA00044	11-06-18
Utah	NELAP	8	CA00044	02-28-18
Virginia	NELAP	3	460278	03-14-18
Washington	State Program	10	C581	05-05-18
West Virginia (DW)	State Program	3	9930C	12-31-18
Wyoming	State Program	8	8TMS-L	01-28-19

Laboratory: TestAmerica Denver

Unless otherwise noted, all analytes for this laboratory were covered under each accreditation/certification below.

Authority	Program	EPA Region	Identification Number	Expiration Date
North Carolina (WW/SW)	State Program	4	358	12-31-18

The following analytes are included in this report, but accreditation/certification is not offered by the governing authority:

Analysis Method	Prep Method	Matrix	Analyte
8321A	3535	Water	HFPO-DA

* Accreditation/Certification renewal pending - accreditation/certification considered valid.

Method Summary

Client: Chemours Company FC, LLC The
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 320-35355-1

Method	Method Description	Protocol	Laboratory
8321A	HFPO-DA	SW846	TAL DEN

Protocol References:

SW846 = "Test Methods For Evaluating Solid Waste, Physical/Chemical Methods", Third Edition, November 1986 And Its Updates.

Laboratory References:

TAL DEN = TestAmerica Denver, 4955 Yarrow Street, Arvada, CO 80002, TEL (303)736-0100

Sample Summary

Client: Chemours Company FC, LLC The
Project/Site: FAY-2018 Residential Sampling

TestAmerica Job ID: 320-35355-1

Lab Sample ID	Client Sample ID	Matrix	Collected	Received
320-35355-1	FAY-D-7154BUTLE-W1-1-012418	Water	01/24/18 15:19	01/25/18 09:30
320-35355-2	FAY-D-384JERML-W1-1-012418	Water	01/24/18 13:35	01/25/18 09:30
320-35355-3	FAY-D-2528INDIA-W1-2-012418	Water	01/24/18 08:56	01/25/18 09:30
320-35355-4	FAY-D-661PINEA-W1-1-012418	Water	01/24/18 09:40	01/25/18 09:30
320-35355-5	FAY-D-701PINEA-W1-1-012418	Water	01/24/18 09:58	01/25/18 09:30
320-35355-6	FAY-D-701PINEA-W1-2-012418	Water	01/24/18 10:05	01/25/18 09:30
320-35355-7	FAY-D-663PINEA-W1-1-012418	Water	01/24/18 10:29	01/25/18 09:30
320-35355-8	FAY-D-7166BUTLE-W2-1-012418	Water	01/24/18 15:43	01/25/18 09:30
320-35355-9	FAY-D-7166BUTLE-W1-1-012418	Water	01/24/18 15:35	01/25/18 09:30
320-35355-10	FAY-D-7166BUTLE-W1-1-012418D	Water	01/24/18 15:35	01/25/18 09:30

TestAmerica Sacramento

LCMS MANUAL INTEGRATION SUMMARY

Lab Name: TestAmerica DenverJob No.: 320-35355-1SDG No.: Instrument ID: LC_LCMS7Analysis Batch Number: 390728Lab Sample ID: STD001 280-390728/3 ICClient Sample ID: Date Analyzed: 10/10/17 09:35Lab File ID: hfpo717J10026.dGC Column: Synergi Hydro ID:

COMPOUND NAME	RETENTION TIME	MANUAL INTEGRATION		
		REASON	ANALYST	DATE
HFPO-DA	0.89	Baseline	meyera	10/10/17 11:50

LCMS MANUAL INTEGRATION SUMMARY

Lab Name: TestAmerica Denver	Job No.: 320-35355-1													
SDG No.:														
Instrument ID: LC_LCMS7	Analysis Batch Number: 403246													
Lab Sample ID: LLCS 280-403139/4-A	Client Sample ID:													
Date Analyzed: 01/30/18 09:25	Lab File ID: hfpo718A30007.d GC Column: Synergi Hydro ID:													
<table border="1"> <thead> <tr> <th rowspan="2">COMPOUND NAME</th> <th rowspan="2">RETENTION TIME</th> <th colspan="3">MANUAL INTEGRATION</th> </tr> <tr> <th>REASON</th> <th>ANALYST</th> <th>DATE</th> </tr> </thead> <tbody> <tr> <td>HFPO-DA</td> <td>0.92</td> <td>Baseline</td> <td>meyera</td> <td>01/30/18 14:47</td> </tr> </tbody> </table>		COMPOUND NAME	RETENTION TIME	MANUAL INTEGRATION			REASON	ANALYST	DATE	HFPO-DA	0.92	Baseline	meyera	01/30/18 14:47
COMPOUND NAME	RETENTION TIME			MANUAL INTEGRATION										
		REASON	ANALYST	DATE										
HFPO-DA	0.92	Baseline	meyera	01/30/18 14:47										
Lab Sample ID: 320-35355-5	Client Sample ID: FAY-D-701PINEA-W1-1-012418													
Date Analyzed: 01/30/18 09:45	Lab File ID: hfpo718A30013.d GC Column: Synergi Hydro ID:													
<table border="1"> <thead> <tr> <th rowspan="2">COMPOUND NAME</th> <th rowspan="2">RETENTION TIME</th> <th colspan="3">MANUAL INTEGRATION</th> </tr> <tr> <th>REASON</th> <th>ANALYST</th> <th>DATE</th> </tr> </thead> <tbody> <tr> <td>HFPO-DA</td> <td>0.97</td> <td>Assign Peak</td> <td>meyera</td> <td>01/30/18 14:48</td> </tr> </tbody> </table>		COMPOUND NAME	RETENTION TIME	MANUAL INTEGRATION			REASON	ANALYST	DATE	HFPO-DA	0.97	Assign Peak	meyera	01/30/18 14:48
COMPOUND NAME	RETENTION TIME			MANUAL INTEGRATION										
		REASON	ANALYST	DATE										
HFPO-DA	0.97	Assign Peak	meyera	01/30/18 14:48										
Lab Sample ID: 320-35355-7	Client Sample ID: FAY-D-663PINEA-W1-1-012418													
Date Analyzed: 01/30/18 09:51	Lab File ID: hfpo718A30015.d GC Column: Synergi Hydro ID:													
<table border="1"> <thead> <tr> <th rowspan="2">COMPOUND NAME</th> <th rowspan="2">RETENTION TIME</th> <th colspan="3">MANUAL INTEGRATION</th> </tr> <tr> <th>REASON</th> <th>ANALYST</th> <th>DATE</th> </tr> </thead> <tbody> <tr> <td>HFPO-DA</td> <td>0.93</td> <td>Baseline</td> <td>meyera</td> <td>01/30/18 14:48</td> </tr> </tbody> </table>		COMPOUND NAME	RETENTION TIME	MANUAL INTEGRATION			REASON	ANALYST	DATE	HFPO-DA	0.93	Baseline	meyera	01/30/18 14:48
COMPOUND NAME	RETENTION TIME			MANUAL INTEGRATION										
		REASON	ANALYST	DATE										
HFPO-DA	0.93	Baseline	meyera	01/30/18 14:48										
Lab Sample ID: 320-35355-9	Client Sample ID: FAY-D-7166BUTLE-W1-1-012418													
Date Analyzed: 01/30/18 09:58	Lab File ID: hfpo718A30017.d GC Column: Synergi Hydro ID:													
<table border="1"> <thead> <tr> <th rowspan="2">COMPOUND NAME</th> <th rowspan="2">RETENTION TIME</th> <th colspan="3">MANUAL INTEGRATION</th> </tr> <tr> <th>REASON</th> <th>ANALYST</th> <th>DATE</th> </tr> </thead> <tbody> <tr> <td>HFPO-DA</td> <td>0.92</td> <td>Baseline</td> <td>meyera</td> <td>02/05/18 13:15</td> </tr> </tbody> </table>		COMPOUND NAME	RETENTION TIME	MANUAL INTEGRATION			REASON	ANALYST	DATE	HFPO-DA	0.92	Baseline	meyera	02/05/18 13:15
COMPOUND NAME	RETENTION TIME			MANUAL INTEGRATION										
		REASON	ANALYST	DATE										
HFPO-DA	0.92	Baseline	meyera	02/05/18 13:15										
Lab Sample ID: 320-35355-9 DU	Client Sample ID: FAY-D-7166BUTLE-W1-1-012418 DU													
Date Analyzed: 01/30/18 10:01	Lab File ID: hfpo718A30018.d GC Column: Synergi Hydro ID:													
<table border="1"> <thead> <tr> <th rowspan="2">COMPOUND NAME</th> <th rowspan="2">RETENTION TIME</th> <th colspan="3">MANUAL INTEGRATION</th> </tr> <tr> <th>REASON</th> <th>ANALYST</th> <th>DATE</th> </tr> </thead> <tbody> <tr> <td>HFPO-DA</td> <td>0.92</td> <td>Baseline</td> <td>meyera</td> <td>02/05/18 13:15</td> </tr> </tbody> </table>		COMPOUND NAME	RETENTION TIME	MANUAL INTEGRATION			REASON	ANALYST	DATE	HFPO-DA	0.92	Baseline	meyera	02/05/18 13:15
COMPOUND NAME	RETENTION TIME			MANUAL INTEGRATION										
		REASON	ANALYST	DATE										
HFPO-DA	0.92	Baseline	meyera	02/05/18 13:15										

REAGENT TRACEABILITY SUMMARY

Lab Name: TestAmerica Denver

Job No.: 320-35355-1

SDG No.:

Reagent ID	Exp Date	Prep Date	Dilutant Used	Reagent Final Volume	Parent Reagent		Analyte	Concentration
					Reagent ID	Volume Added		
HFPO I.S._00007	12/12/18	12/12/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA_00007	1 mL	13C3 HFPO-DA	0.5 ug/mL
							13C3 HFPO-DA (IS)	0.5 ug/mL
.13C3 HFPO-DA_00007	08/17/20	Wellington Laboratories, Lot M3HFPOADA0817		(Purchased Reagent)		13C3 HFPO-DA	50 ug/mL	
HFPO Spike_00004	10/30/18	10/30/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	HFPO-DA_00004	1 mL	HFPO-DA	0.5 ug/mL
							HFPO-DA	50 ug/mL
HFPO_CAL-0_00031	10/24/17	10/10/17	PFC Dill_Solvent, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	13C3 HFPO-DA	0.5 ug/mL
.13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616		(Purchased Reagent)		13C3 HFPO-DA	50 ug/mL	
HFPO_CAL-1_00030	09/28/17	09/14/17	80:20 Methanol : H2O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L
							13C3 HFPO-DA (IS)	10 ug/L
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	HFPO-DA	0.25 ug/L
.13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616		(Purchased Reagent)		13C3 HFPO-DA	0.5 ug/mL	
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	HFPO-DA_00003	1 mL	13C3 HFPO-DA (IS)	0.5 ug/mL
.HFPO-DA_00003	12/16/18	Wellington Laboratories, Lot HFPOADA0213		(Purchased Reagent)		HFPO-DA	HFPO-DA	50 ug/mL
HFPO_CAL-1_00031	10/24/17	10/10/17	80:20 Methanol : H2O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	0.5 ug/mL
							13C3 HFPO-DA (IS)	0.25 ug/L
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	HFPO-DA	0.5 ug/mL
.13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616		(Purchased Reagent)		13C3 HFPO-DA	50 ug/mL	
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	HFPO-DA_00003	1 mL	13C3 HFPO-DA (IS)	0.5 ug/mL
.HFPO-DA_00003	12/16/18	Wellington Laboratories, Lot HFPOADA0213		(Purchased Reagent)		HFPO-DA	HFPO-DA	50 ug/mL
HFPO_CAL-2_00031	09/28/17	09/14/17	80:20 Methanol : H2O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L
							13C3 HFPO-DA (IS)	10 ug/L
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	HFPO-DA	0.5 ug/mL
.13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616		(Purchased Reagent)		13C3 HFPO-DA (IS)	0.5 ug/mL	

REAGENT TRACEABILITY SUMMARY

Lab Name: TestAmerica Denver

Job No.: 320-35355-1

SDG No.:

Reagent ID	Exp Date	Prep Date	Dilutant Used	Reagent Final Volume	Parent Reagent		Analyte	Concentration	
					Reagent ID	Volume Added			
..13C3 HFPO-DA_00004	08/28/18		Wellington Laboratories, Lot M3HFPOADA0616		(Purchased Reagent)		13C3 HFPO-DA	50 ug/mL	
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	HFPO-DA_00003	1 mL	HFPO-DA	0.5 ug/mL	
..HFPO-DA_00003	12/16/18		Wellington Laboratories, Lot HFPOADA0213		(Purchased Reagent)		HFPO-DA	50 ug/mL	
HFPO_CAL-2_00032	10/24/17	10/10/17	80:20 Methanol : H ₂ O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L	
					HFPO Spike_00003	1 uL	13C3 HFPO-DA (IS)	10 ug/L	
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	HFPO-DA	0.5 ug/mL	
							13C3 HFPO-DA (IS)	0.5 ug/mL	
..13C3 HFPO-DA_00004	08/28/18		Wellington Laboratories, Lot M3HFPOADA0616		(Purchased Reagent)		13C3 HFPO-DA	50 ug/mL	
							13C3 HFPO-DA (IS)	50 ug/mL	
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	HFPO-DA_00003	1 mL	HFPO-DA	0.5 ug/mL	
							HFPO-DA	50 ug/mL	
..HFPO-DA_00003	12/16/18		Wellington Laboratories, Lot HFPOADA0213		(Purchased Reagent)		HFPO-DA	50 ug/mL	
							HFPO-DA	50 ug/mL	
HFPO_CAL-3_00030	09/28/17	09/14/17	80:20 Methanol : H ₂ O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L	
					HFPO Spike_00003	2 uL	13C3 HFPO-DA (IS)	10 ug/L	
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	HFPO-DA	0.5 ug/mL	
							13C3 HFPO-DA (IS)	0.5 ug/mL	
..13C3 HFPO-DA_00004	08/28/18		Wellington Laboratories, Lot M3HFPOADA0616		(Purchased Reagent)		13C3 HFPO-DA	50 ug/mL	
							13C3 HFPO-DA (IS)	50 ug/mL	
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	HFPO-DA_00003	1 mL	HFPO-DA	0.5 ug/mL	
							HFPO-DA	50 ug/mL	
..HFPO-DA_00003	12/16/18		Wellington Laboratories, Lot HFPOADA0213		(Purchased Reagent)		HFPO-DA	50 ug/mL	
							HFPO-DA	50 ug/mL	
HFPO_CAL-3_00031	10/24/17	10/10/17	80:20 Methanol : H ₂ O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L	
					HFPO Spike_00003	2 uL	13C3 HFPO-DA (IS)	10 ug/L	
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	HFPO-DA	0.5 ug/mL	
							13C3 HFPO-DA (IS)	0.5 ug/mL	
..13C3 HFPO-DA_00004	08/28/18		Wellington Laboratories, Lot M3HFPOADA0616		(Purchased Reagent)		13C3 HFPO-DA	50 ug/mL	
							13C3 HFPO-DA (IS)	50 ug/mL	
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	HFPO-DA_00003	1 mL	HFPO-DA	0.5 ug/mL	
							HFPO-DA	50 ug/mL	
..HFPO-DA_00003	12/16/18		Wellington Laboratories, Lot HFPOADA0213		(Purchased Reagent)		HFPO-DA	50 ug/mL	
							HFPO-DA	50 ug/mL	
HFPO_CAL-4_00030	09/28/17	09/14/17	80:20 Methanol : H ₂ O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L	
							13C3 HFPO-DA (IS)	10 ug/L	

REAGENT TRACEABILITY SUMMARY

Lab Name: TestAmerica Denver

Job No.: 320-35355-1

SDG No.:

Reagent ID	Exp Date	Prep Date	Dilutant Used	Reagent Final Volume	Parent Reagent		Analyte	Concentration
					Reagent ID	Volume Added		
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	HFPO Spike_00003	4 uL	HFPO-DA	2 ug/L
					13C3 HFPO-DA_00004	1 mL	13C3 HFPO-DA	0.5 ug/mL
..13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616			(Purchased Reagent)		13C3 HFPO-DA (IS)	0.5 ug/mL
							13C3 HFPO-DA	50 ug/mL
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	HFPO-DA_00003	1 mL	HFPO-DA (IS)	50 ug/mL
..HFPO-DA_00003	12/16/18	Wellington Laboratories, Lot HFPOADA0213			(Purchased Reagent)		HFPO-DA	0.5 ug/mL
							HFPO-DA	50 ug/mL
HFPO_CAL-4_00031	10/24/17	10/10/17	80:20 Methanol : H ₂ O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L
							13C3 HFPO-DA (IS)	10 ug/L
					HFPO Spike_00003	4 uL	HFPO-DA	2 ug/L
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	13C3 HFPO-DA	0.5 ug/mL
..13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616			(Purchased Reagent)		13C3 HFPO-DA (IS)	0.5 ug/mL
							13C3 HFPO-DA	50 ug/mL
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	HFPO-DA_00003	1 mL	HFPO-DA	0.5 ug/mL
..HFPO-DA_00003	12/16/18	Wellington Laboratories, Lot HFPOADA0213			(Purchased Reagent)		HFPO-DA	0.5 ug/mL
							HFPO-DA	50 ug/mL
HFPO_CAL-5_00067	09/28/17	09/14/17	80:20 Methanol : H ₂ O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L
							13C3 HFPO-DA (IS)	10 ug/L
					HFPO Spike_00003	10 uL	HFPO-DA	5 ug/L
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	13C3 HFPO-DA	0.5 ug/mL
..13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616			(Purchased Reagent)		13C3 HFPO-DA (IS)	0.5 ug/mL
							13C3 HFPO-DA	50 ug/mL
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	HFPO-DA_00003	1 mL	HFPO-DA	0.5 ug/mL
..HFPO-DA_00003	12/16/18	Wellington Laboratories, Lot HFPOADA0213			(Purchased Reagent)		HFPO-DA	0.5 ug/mL
							HFPO-DA	50 ug/mL
HFPO_CAL-5_00070	10/24/17	10/10/17	80:20 Methanol : H ₂ O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L
							13C3 HFPO-DA (IS)	10 ug/L
					HFPO Spike_00003	10 uL	HFPO-DA	5 ug/L
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	13C3 HFPO-DA	0.5 ug/mL
..13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616			(Purchased Reagent)		13C3 HFPO-DA (IS)	0.5 ug/mL
							13C3 HFPO-DA	50 ug/mL
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	HFPO-DA_00003	1 mL	HFPO-DA	0.5 ug/mL
..HFPO-DA_00003	12/16/18	Wellington Laboratories, Lot HFPOADA0213			(Purchased Reagent)		HFPO-DA	0.5 ug/mL
							HFPO-DA	50 ug/mL

REAGENT TRACEABILITY SUMMARY

Lab Name: TestAmerica Denver

Job No.: 320-35355-1

SDG No.:

Reagent ID	Exp Date	Prep Date	Dilutant Used	Reagent Final Volume	Parent Reagent		Analyte	Concentration
					Reagent ID	Volume Added		
HFPO_CAL-5_00079	02/09/18	01/26/18	80:20 Methanol : H ₂ O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L
					HFPO Spike_00004	10 uL	HFPO-DA	5 ug/L
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	13C3 HFPO-DA	0.5 ug/mL
..13C3 HFPO-DA_00004	08/28/18		Wellington Laboratories, Lot M3HFPOADA0616		(Purchased Reagent)		13C3 HFPO-DA	50 ug/mL
.HFPO Spike_00004	10/30/18	10/30/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	HFPO-DA_00004	1 mL	HFPO-DA	0.5 ug/mL
..HFPO-DA_00004	07/13/20		Wellington Laboratories, Lot HFPOADA0717		(Purchased Reagent)		HFPO-DA	50 ug/mL
HFPO_CAL-6_00067	09/28/17	09/14/17	80:20 Methanol : H ₂ O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L
					HFPO Spike_00003	20 uL	13C3 HFPO-DA (IS)	10 ug/L
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	HFPO-DA	0.5 ug/mL
..13C3 HFPO-DA_00004	08/28/18		Wellington Laboratories, Lot M3HFPOADA0616		(Purchased Reagent)		13C3 HFPO-DA	50 ug/mL
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	HFPO-DA_00003	1 mL	HFPO-DA	0.5 ug/mL
..HFPO-DA_00003	12/16/18		Wellington Laboratories, Lot HFPOADA0213		(Purchased Reagent)		HFPO-DA	50 ug/mL
HFPO_CAL-6_00070	10/24/17	10/10/17	80:20 Methanol : H ₂ O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L
					HFPO Spike_00003	20 uL	13C3 HFPO-DA (IS)	10 ug/L
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	HFPO-DA	0.5 ug/mL
..13C3 HFPO-DA_00004	08/28/18		Wellington Laboratories, Lot M3HFPOADA0616		(Purchased Reagent)		13C3 HFPO-DA	50 ug/mL
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	HFPO-DA_00003	1 mL	HFPO-DA	0.5 ug/mL
..HFPO-DA_00003	12/16/18		Wellington Laboratories, Lot HFPOADA0213		(Purchased Reagent)		HFPO-DA	50 ug/mL
HFPO_CAL-6_00079	02/09/18	01/26/18	80:20 Methanol : H ₂ O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L
					HFPO Spike_00004	20 uL	HFPO-DA	10 ug/L
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	13C3 HFPO-DA	0.5 ug/mL
..13C3 HFPO-DA_00004	08/28/18		Wellington Laboratories, Lot M3HFPOADA0616		(Purchased Reagent)		13C3 HFPO-DA	50 ug/mL
.HFPO Spike_00004	10/30/18	10/30/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	HFPO-DA_00004	1 mL	HFPO-DA	0.5 ug/mL
..HFPO-DA_00004	07/13/20		Wellington Laboratories, Lot HFPOADA0717		(Purchased Reagent)		HFPO-DA	50 ug/mL
HFPO_CAL-7_00030	09/28/17	09/14/17	80:20 Methanol : H ₂ O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L

REAGENT TRACEABILITY SUMMARY

Lab Name: TestAmerica Denver

Job No.: 320-35355-1

SDG No.:

Reagent ID	Exp Date	Prep Date	Dilutant Used	Reagent Final Volume	Parent Reagent		Analyte	Concentration			
					Reagent ID	Volume Added					
.HFPO_I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	13C3 HFPO-DA (IS)	10 ug/L			
							HFPO-DA	25 ug/L			
..13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616			(Purchased Reagent)		13C3 HFPO-DA	0.5 ug/mL			
							13C3 HFPO-DA (IS)	0.5 ug/mL			
.HFPO_Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	HFPO-DA_00003	1 mL	HFPO-DA	50 ug/mL			
..HFPO-DA_00003	12/16/18	Wellington Laboratories, Lot HFPOADA0213			(Purchased Reagent)		HFPO-DA	50 ug/mL			
HFPO_CAL-7_00031	10/24/17	10/10/17	80:20 Methanol : H2O, Lot 00016	1 mL	HFPO_I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L			
							13C3 HFPO-DA (IS)	10 ug/L			
.HFPO_I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	HFPO-DA	25 ug/L			
							13C3 HFPO-DA	0.5 ug/mL			
..13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616			(Purchased Reagent)		13C3 HFPO-DA (IS)	0.5 ug/mL			
							13C3 HFPO-DA	50 ug/mL			
.HFPO_Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	HFPO-DA_00003	1 mL	HFPO-DA	0.5 ug/mL			
..HFPO-DA_00003	12/16/18	Wellington Laboratories, Lot HFPOADA0213			(Purchased Reagent)		HFPO-DA	50 ug/mL			
HFPO_CAL-8_00030	09/28/17	09/14/17	80:20 Methanol : H2O, Lot 00016	1 mL	HFPO_I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L			
							13C3 HFPO-DA (IS)	10 ug/L			
.HFPO_I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	HFPO-DA	50 ug/L			
							13C3 HFPO-DA	0.5 ug/mL			
..13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616			(Purchased Reagent)		13C3 HFPO-DA (IS)	0.5 ug/mL			
							13C3 HFPO-DA	50 ug/mL			
.HFPO_Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	HFPO-DA_00003	1 mL	HFPO-DA	0.5 ug/mL			
..HFPO-DA_00003	12/16/18	Wellington Laboratories, Lot HFPOADA0213			(Purchased Reagent)		HFPO-DA	50 ug/mL			
HFPO_CAL-8_00031	10/24/17	10/10/17	80:20 Methanol : H2O, Lot 00016	1 mL	HFPO_I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L			
							13C3 HFPO-DA (IS)	10 ug/L			
.HFPO_I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	HFPO-DA	50 ug/L			
							13C3 HFPO-DA	0.5 ug/mL			
..13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616			(Purchased Reagent)		13C3 HFPO-DA (IS)	0.5 ug/mL			
							13C3 HFPO-DA	50 ug/mL			
.HFPO_Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS_MeOH_00110	100 mL	HFPO-DA_00003	1 mL	HFPO-DA	0.5 ug/mL			

REAGENT TRACEABILITY SUMMARY

Lab Name: TestAmerica Denver

Job No.: 320-35355-1

SDG No.:

Reagent ID	Exp Date	Prep Date	Dilutant Used	Reagent Final Volume	Parent Reagent		Analyte	Concentration
					Reagent ID	Volume Added		
..HFPO-DA_00003	12/16/18	Wellington Laboratories, Lot HFPODA0213		(Purchased Reagent)		HFPO-DA	50 ug/mL	
HFPO_ICV_00031	09/28/17	09/14/17	80:20 Methanol : H ₂ O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L
					HFPO Spike_00003	4 uL		
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	13C3 HFPO-DA	0.5 ug/mL
..13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616		(Purchased Reagent)		13C3 HFPO-DA	50 ug/mL	
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS MeOH_00110	100 mL	HFPO-DA_00003	1 mL	HFPO-DA	0.5 ug/mL
..HFPO-DA_00003	12/16/18	Wellington Laboratories, Lot HFPODA0213		(Purchased Reagent)		HFPO-DA	50 ug/mL	
HFPO_ICV_00032	10/24/17	10/10/17	80:20 Methanol : H ₂ O, Lot 00016	1 mL	HFPO I.S._00004	20 uL	13C3 HFPO-DA	10 ug/L
					HFPO Spike_00003	4 uL		
.HFPO I.S._00004	08/28/18	08/28/17	LCMS Grade MeOH, Lot LCMS MeOH_00110	100 mL	13C3 HFPO-DA_00004	1 mL	13C3 HFPO-DA	0.5 ug/mL
..13C3 HFPO-DA_00004	08/28/18	Wellington Laboratories, Lot M3HFPOADA0616		(Purchased Reagent)		13C3 HFPO-DA	50 ug/mL	
.HFPO Spike_00003	01/11/18	01/10/17	LCMS Grade MeOH, Lot LCMS MeOH_00110	100 mL	HFPO-DA_00003	1 mL	HFPO-DA	0.5 ug/mL
..HFPO-DA_00003	12/16/18	Wellington Laboratories, Lot HFPODA0213		(Purchased Reagent)		HFPO-DA	50 ug/mL	

Reagent

13C3 HFPO-DA_00004



WELLINGTON LABORATORIES

CERTIFICATE OF ANALYSIS DOCUMENTATION


PRODUCT CODE: M3HFPO-DA

LOT NUMBER: M3HFPODA0616

COMPOUND: 2,3,3,3-Tetrafluoro-2-(1,1,2,2,3,3,3-heptafluoropropoxy)-¹³C₃-propanoic acid

STRUCTURE:

CAS #: Not available

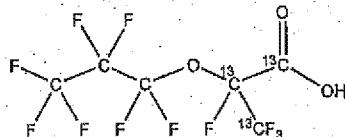

MOLECULAR FORMULA: ¹³C₃¹²C₃HF₁₁O₃

MOLECULAR WEIGHT: 333.03

CONCENTRATION: 50 ± 2.5 µg/ml

SOLVENT(S): Methanol

CHEMICAL PURITY: >98%

ISOTOPIC PURITY: >99% ¹³C

LAST TESTED: (mm/dd/yyyy) 06/25/2016

(¹³C₃)

EXPIRY DATE: (mm/dd/yyyy) 06/25/2019

RECOMMENDED STORAGE: Store ampoule in a cool, dark place

DOCUMENTATION/ DATA ATTACHED:

Figure 1: LC/MS Data (TIC and Mass Spectrum)

Figure 2: LC/MS/MS Data (Selected MRM Transitions)

ADDITIONAL INFORMATION:

- See page 2 for further details.
- Contains ~ 1.5% of two constitutional isomers.
- Product is commercially known as GenX.

FOR LABORATORY USE ONLY: NOT FOR HUMAN OR DRUG USE

Certified By:

B.G. Chittim

Date: 06/29/2016

(mm/dd/yyyy)

Wellington Laboratories Inc., 345 Southgate Dr. Guelph ON N1G 3M5 CANADA
519-822-2436 • Fax: 519-822-2849 • info@well-labs.com

INTENDED USE:

The products prepared by Wellington Laboratories Inc. are for laboratory use only. This certified reference material (CRM) was designed to be used as a standard for the identification and/or quantification of the specific chemical compound it contains.

HAZARDS:

This product should only be used by qualified personnel familiar with its potential hazards and trained in the handling of hazardous chemicals. Due care should be exercised to prevent unnecessary human contact or ingestion. All procedures should be carried out in a well-functioning fume hood and suitable gloves, eye protection, and clothing should be worn at all times. Waste should be disposed of according to national and regional regulations. Safety Data Sheets (SDSs) are available upon request.

SYNTHESIS / CHARACTERIZATION:

Where possible, all of our products are synthesized using single-product unambiguous routes. They are then characterized, and their structures and purities confirmed, using a combination of the most relevant techniques, such as NMR, GC/MS, LC/MS/MS, SFC/UV/MS/MS, x-ray crystallography, and melting point. Isotopic purities of mass-labelled compounds are also confirmed using HRGC/HRMS and/or LC/MS/MS.

HOMOGENEITY:

Prior to solution preparation, crystalline material is tested for homogeneity using a variety of techniques (as stated above) and its solubility in a given diluent is taken into consideration. Duplicate solutions of a new product are prepared from the same crystalline lot and, after the addition of an appropriate internal standard, they are compared by GC/MS, LC/MS/MS and/or SFC/UV/MS/MS. The relative response factors of the analyte of interest in each solution are required to be <5% RSD. New solution lots of existing products are compared to older lots in the same manner, which further confirms the homogeneity of the crystalline material as well as the stability and homogeneity of the solutions in the storage containers.

UNCERTAINTY:

The maximum combined relative standard uncertainty of our reference standard solutions is calculated using the following equation:

The combined relative standard uncertainty, $u_c(y)$, of a value y and the uncertainty of the independent parameters x_1, x_2, \dots, x_n on which it depends is:

$$u_c(y(x_1, x_2, \dots, x_n)) = \sqrt{\sum_{i=1}^n u(y, x_i)^2}$$

where x is expressed as a relative standard uncertainty of the individual parameter.

The individual uncertainties taken into account include those associated with weights (calibration of the balance) and volumes (calibration of the volumetric glassware). An expanded maximum combined percent relative uncertainty of ±5% (calculated with a coverage factor of 2 and a level of confidence of 95%) is stated on the Certificate of Analysis for all of our products.

TRACEABILITY:

All reference standard solutions are traceable to specific crystalline lots. The microbalances used for solution preparation are regularly tested by an external ISO/IEC 17025 accredited calibration company. In addition, their calibration is verified prior to each weighing using NIST and/or NRC traceable external weights. All volumetric glassware used is of Class A tolerance and has been tested according to the appropriate ASTM procedures, which are ultimately traceable to NIST. For certain products, traceability to international interlaboratory studies has also been established.

EXPIRY DATE / PERIOD OF VALIDITY:

Ongoing stability studies of this product have demonstrated stability in its composition and concentration, until the specified expiry date, in the unopened ampoule. Monitoring for any degradation or change in concentration of the listed analyte(s) is performed on a routine basis.

LIMITED WARRANTY:

At the time of shipment, all products are warranted to be free of defects in material and workmanship and to conform to the stated technical and purity specifications.



QUALITY MANAGEMENT:

This product was produced using a Quality Management System registered to the latest versions of ISO 9001 by SAI Global, ISO/IEC 17025 by the Canadian Association for Laboratory Accreditation Inc. (CALA; A 1226), and ISO GUIDE 34 by ANSI-ASQ National Accreditation Board (ANAB; AR-1523).



For additional information or assistance concerning this or any other products from Wellington Laboratories Inc., please visit our website at www.well-labs.com or contact us directly at info@well-labs.com

Figure 1: M3HFPO-DA; LC/MS Data (TIC and Mass Spectrum)

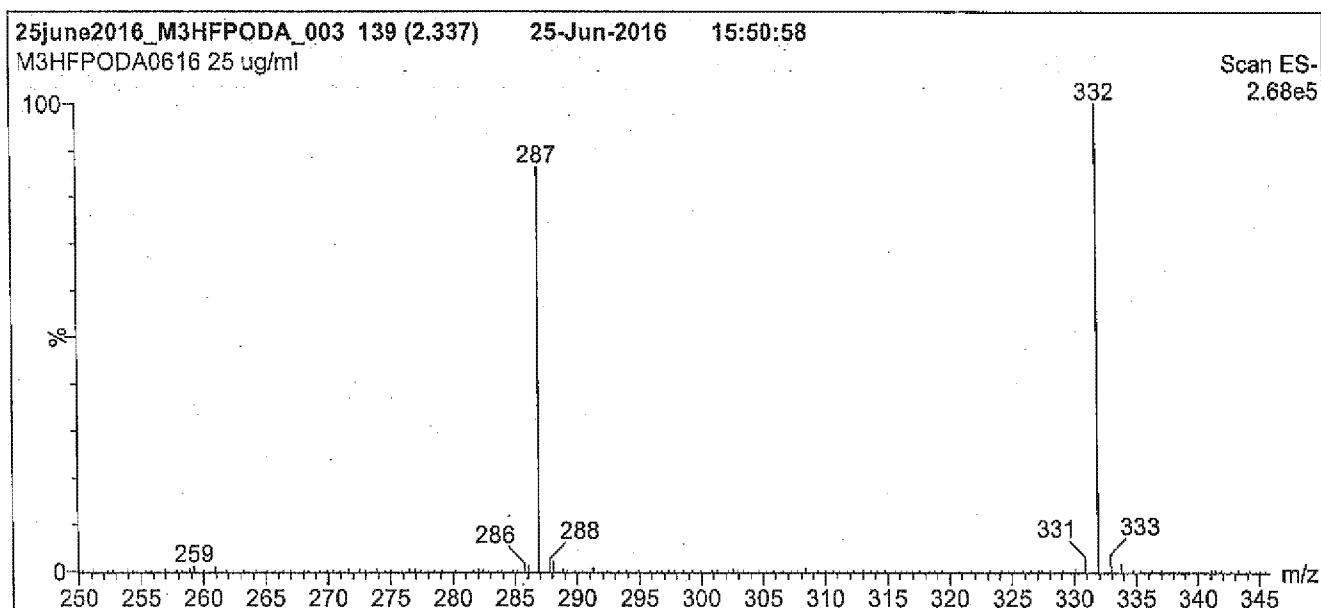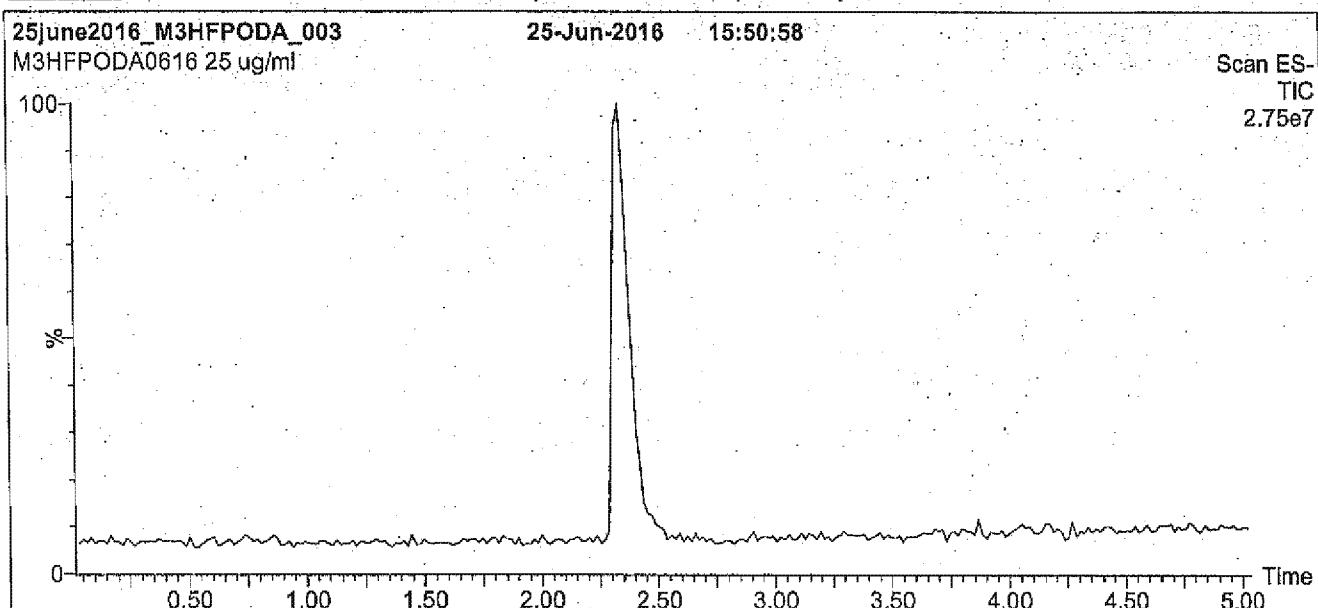

Conditions for Figure 1:

LC: Waters Acquity Ultra Performance LC
MS: Micromass Quattro micro API MS

Chromatographic Conditions

Column: Acquity UPLC BEH Shield RP₁₈,
1.7 μ m, 2.1 x 100 mm

Mobile phase: Gradient
Start: 50% (80:20 MeOH:ACN) / 50% H₂O
(both with 10 mM NH₄OAc buffer)
Ramp to 90% organic over 7 min and hold for 1.5 min
before returning to initial conditions in 0.5 min.
Time: 10 min


Flow: 300 μ l/min

MS Parameters

Experiment: Full Scan (250 - 850 amu)

Source: Electrospray (negative)
Capillary Voltage (kV) = 3.00
Cone Voltage (V) = 9.00
Cone Gas Flow (l/hr) = 50
Desolvation Gas Flow (l/hr) = 750

Figure 2: M3HFPO-DA; LC/MS/MS Data (Selected MRM Transitions)

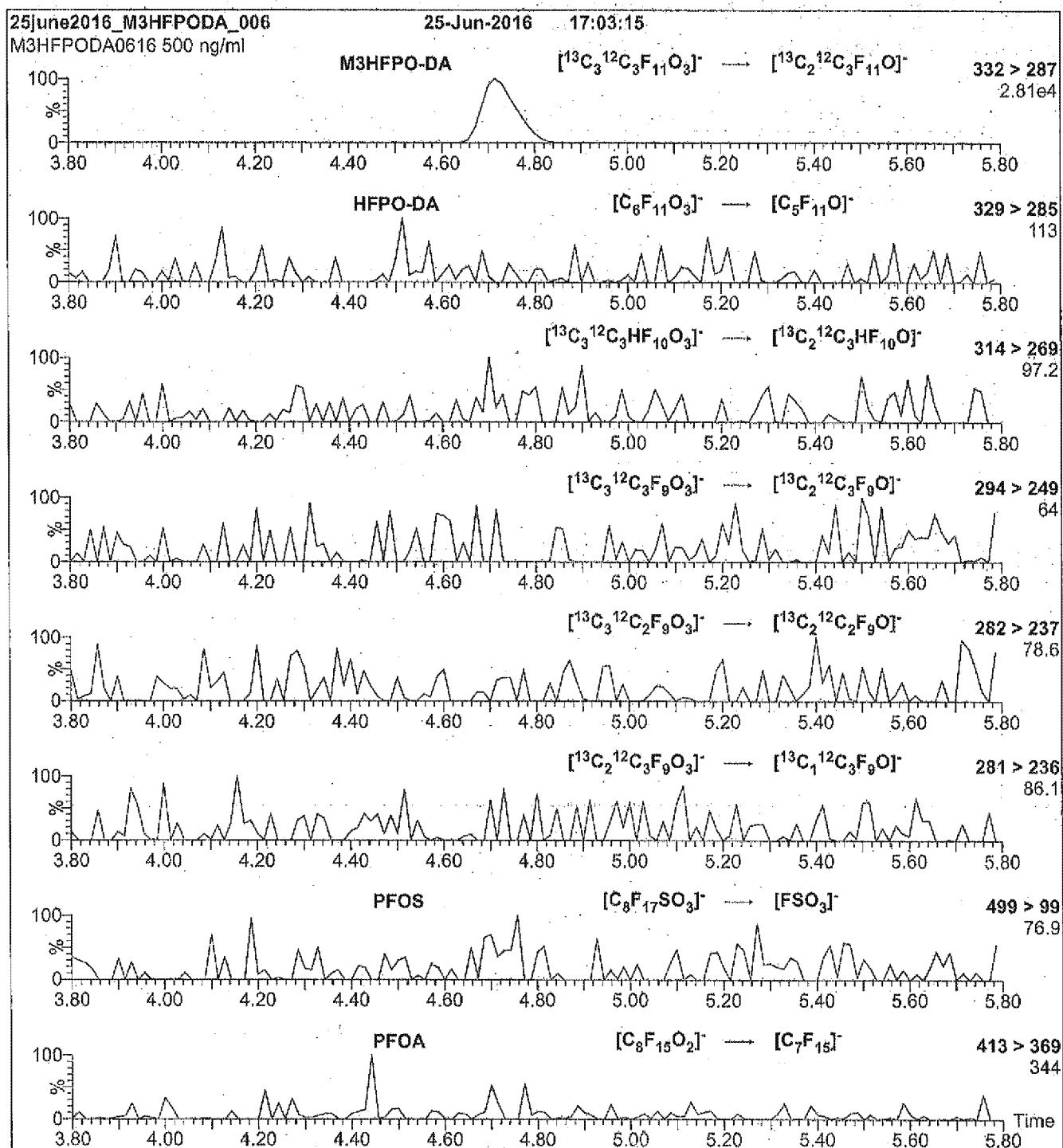

Conditions for Figure 2:

Injection: Direct loop injection
10 μ l (500 ng/ml M3HFPO-DA)

MS Parameters

Collision Gas (mbar) = 3.46e-3
Collision Energy (eV) = 5

Mobile phase: Isocratic 80% (80:20 MeOH:ACN) / 20% H₂O
(both with 10 mM NH₄OAc buffer)

Flow: 300 μ l/min

Reagent

13C3 HFPO-DA_00007



WELLINGTON LABORATORIES

CERTIFICATE OF ANALYSIS DOCUMENTATION

PRODUCT CODE:


M3HFPO-DA

LOT NUMBER: M3HFPODA0817

COMPOUND:

2,3,3,3-Tetrafluoro-2-(1,1,2,2,3,3-heptafluoropropoxy)-¹³C₃-propanoic acid

STRUCTURE:

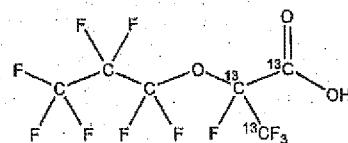

CAS #:

Not available.

MOLECULAR FORMULA:

¹³C₃¹²C₃HF₁₁O₃

MOLECULAR WEIGHT: 333.03

CONCENTRATION:

50 ± 2.5 µg/ml

Methanol

CHEMICAL PURITY:

>98%

SOLVENT(S):

LAST TESTED: (mm/dd/yyyy)

08/17/2017

ISOTOPIC PURITY: >99% ¹³C
(¹³C₃)

EXPIRY DATE: (mm/dd/yyyy)

08/17/2020

RECOMMENDED STORAGE:

Store ampoule in a cool, dark place

DOCUMENTATION/DATA ATTACHED:

Figure 1: LC/MS Data (TIC and Mass Spectrum)

Figure 2: LC/MS/MS Data (Selected MRM Transitions)

ADDITIONAL INFORMATION:

- See page 2 for further details.
- Contains ~ 1.5% of two constitutional isomers.
- Product is commercially known as GenX.

FOR LABORATORY USE ONLY: NOT FOR HUMAN OR DRUG USE

Certified By:

B.G. Chittim, General Manager

Date: 08/25/2017

(mm/dd/yyyy)

INTENDED USE:

The products prepared by Wellington Laboratories Inc. are for laboratory use only. This certified reference material (CRM) was designed to be used as a standard for the identification and/or quantification of the specific chemical compound it contains.

HAZARDS:

This product should only be used by qualified personnel familiar with its potential hazards and trained in the handling of hazardous chemicals. Due care should be exercised to prevent unnecessary human contact or ingestion. All procedures should be carried out in a well-functioning fume hood and suitable gloves, eye protection, and clothing should be worn at all times. Waste should be disposed of according to national and regional regulations. Safety Data Sheets (SDSs) are available upon request.

SYNTHESIS / CHARACTERIZATION:

Where possible, all of our products are synthesized using single-product unambiguous routes. They are then characterized, and their structures and purities confirmed, using a combination of the most relevant techniques, such as NMR, GC/MS, LC/MS/MS, SFC/UW/MS/MS, x-ray crystallography, and melting point. Isotopic purities of mass-labelled compounds are also confirmed using HRGC/HRMS and/or LC/MS/MS.

HOMOGENEITY:

Prior to solution preparation, crystalline material is tested for homogeneity using a variety of techniques (as stated above) and its solubility in a given diluent is taken into consideration. Duplicate solutions of a new product are prepared from the same crystalline lot and, after the addition of an appropriate internal standard, they are compared by GC/MS, LC/MS/MS and/or SFC/UW/MS/MS. The relative response factors of the analyte of interest in each solution are required to be <5% RSD. New solution lots of existing products are compared to older lots in the same manner, which further confirms the homogeneity of the crystalline material as well as the stability and homogeneity of the solutions in the storage containers. In order to maintain the integrity of the assigned value(s), and associated uncertainty, the dilution or injection of a subsample of this product should be performed using calibrated measuring equipment.

UNCERTAINTY:

The maximum combined relative standard uncertainty of our reference standard solutions is calculated using the following equation:

The combined relative standard uncertainty, $u_c(y)$, of a value y and the uncertainty of the independent parameters

x_1, x_2, \dots, x_n on which it depends is:

$$u_c(y(x_1, x_2, \dots, x_n)) = \sqrt{\sum_{i=1}^n u(y, x_i)^2}$$

where x is expressed as a relative standard uncertainty of the individual parameter.

The individual uncertainties taken into account include those associated with weights (calibration of the balance) and volumes (calibration of the volumetric glassware). An expanded maximum combined percent relative uncertainty of ±5% (calculated with a coverage factor of 2 and a level of confidence of 95%) is stated on the Certificate of Analysis for all of our products.

TRACEABILITY:

All reference standard solutions are traceable to specific crystalline lots. The microbalances used for solution preparation are regularly tested by an external ISO/IEC 17025 accredited calibration company. In addition, their calibration is verified prior to each weighing using calibrated NIST and/or NRC traceable external weights. All volumetric glassware used is calibrated, of Class A tolerance, and has been tested according to the appropriate ASTM procedures, which are ultimately traceable to NIST. For certain products, traceability to international interlaboratory studies has also been established.

EXPIRY DATE / PERIOD OF VALIDITY:

Ongoing stability studies of this product have demonstrated stability in its composition and concentration, until the specified expiry date, in the unopened ampoule. Monitoring for any degradation or change in concentration of the listed analyte(s) is performed on a routine basis.

LIMITED WARRANTY:

At the time of shipment, all products are warranted to be free of defects in material and workmanship and to conform to the stated technical and purity specifications.



QUALITY MANAGEMENT:

This product was produced using a Quality Management System registered to the latest versions of ISO 9001 by SAI Global, ISO/IEC 17025 by the Canadian Association for Laboratory Accreditation Inc. (CALA; A 1226), and ISO GUIDE 34 by ANSI-ASQ National Accreditation Board (ANAB; AR-1523).



For additional information or assistance concerning this or any other products from Wellington Laboratories Inc., please visit our website at www.well-labs.com or contact us directly at info@well-labs.com

Figure 1: M3HFPO-DA; LC/MS Data (TIC and Mass Spectrum)

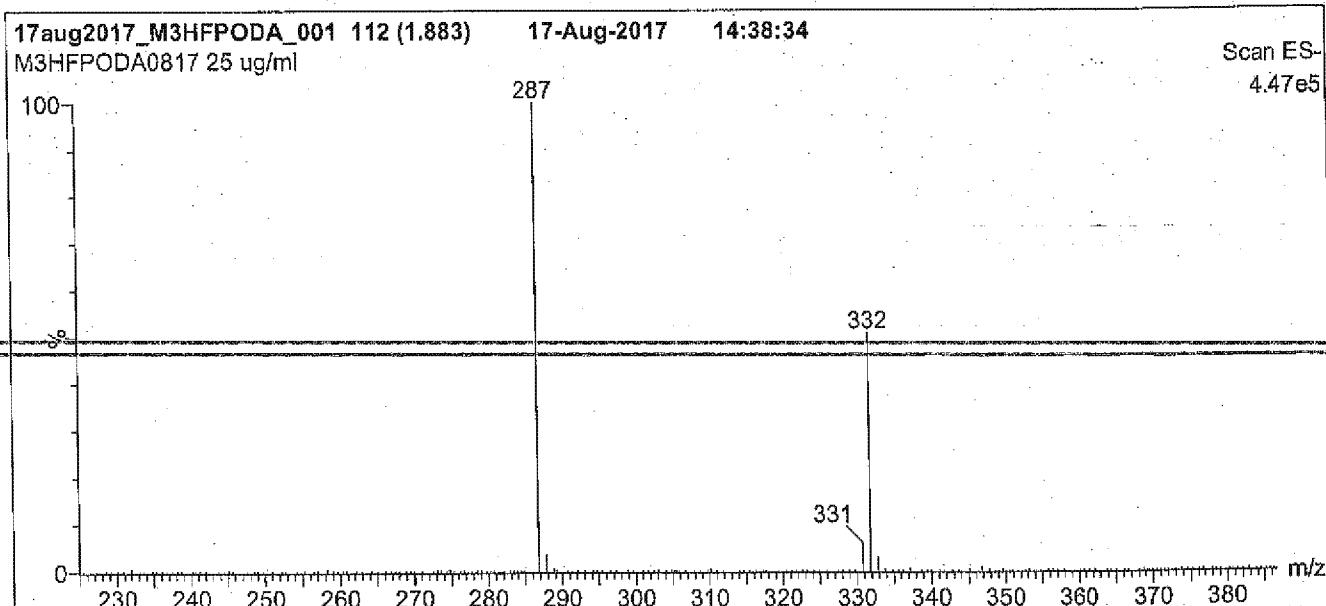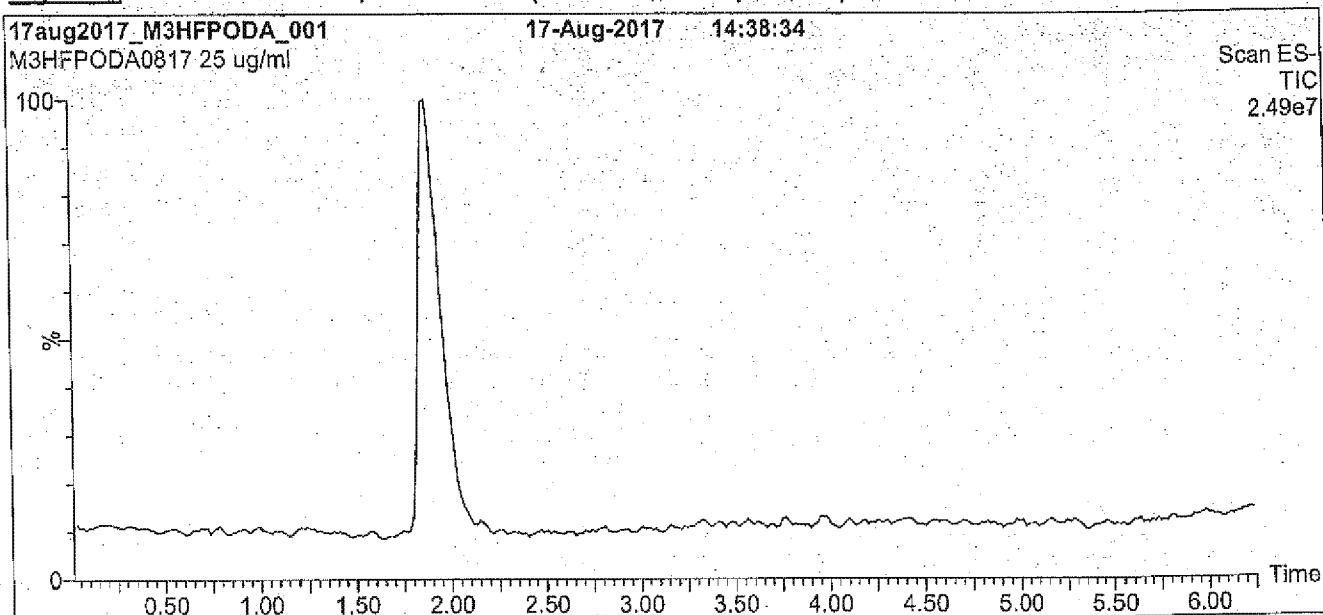

Conditions for Figure 1:

LC: Waters Acquity Ultra Performance LC
MS: Micromass Quattro *micro* API MS

Chromatographic Conditions

Column: Acquity UPLC BEH Shield RP₁₈
 1.7 μm, 2.1 x 100 mm

Mobile phase: Gradient

Start: 55% MeOH / 45% H₂O with 10 mM NH₄OAc buffer
 Ramp to 90% organic over 7.5 min and hold for 1.5 min
 before returning to initial conditions in 0.5 min.

Time: 10 min


Flow: 300 μl/min

MS Parameters

Experiment: Full Scan (225 - 850 amu)

Source: Electrospray (negative)
 Capillary Voltage (kV) = 3.00
 Cone Voltage (V) = 10.00
 Cone Gas Flow (l/hr) = 100
 Desolvation Gas Flow (l/hr) = 750

Figure 2: M3HFPO-DA; LC/MS/MS Data (Selected MRM Transitions)

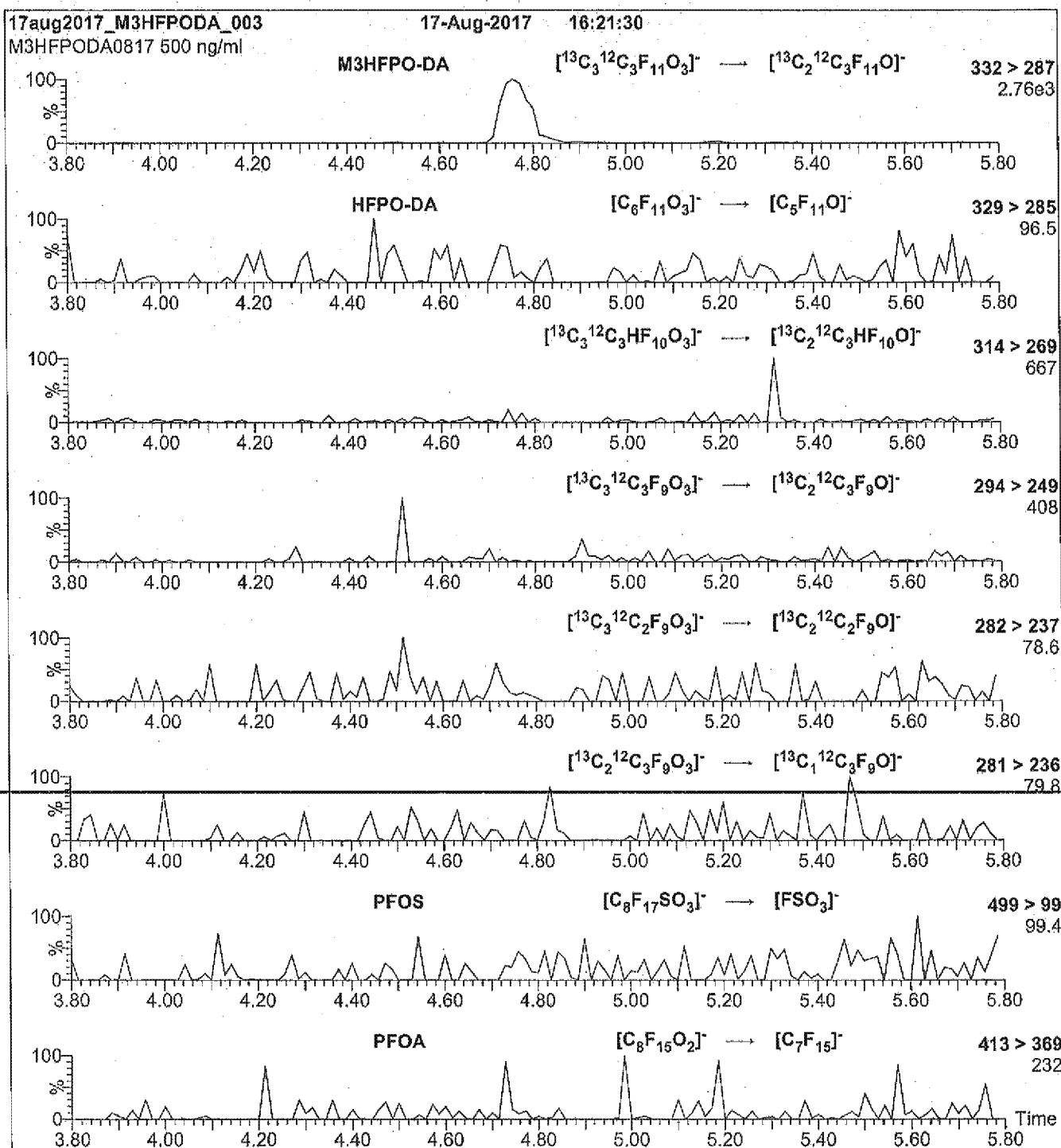

Conditions for Figure 2:

Injection: Direct loop injection
10 μ l (500 ng/ml M3HFPO-DA)

MS Parameters

Collision Gas (mbar) = 3.24e-3
Collision Energy (eV) = 5

Mobile phase: Isocratic 80% MeOH / 20% H₂O wth 10 mM NH₄OAc buffer

Flow: 300 μ l/min

Reagent

HFPO I.S._00004

**Reagent ID: HFPO I.S._00004**

Description:	Internal Standard for HFPO 0.5ug/ml	Expiration Date:	08/28/2018
No. of Bottles:	1	Laboratory:	TestAmerica Denver
Storage Location:	North Analytical	Prepared By:	Meyer, Andrew GC
Reagent Volume:	100.000 mL	Solvent:	LCMS Grade MeOH
Creation Date:	08/28/2017	Solvent Lot#:	LCMS_MeOH_00110
Open Date:			
Container(s):	4700620		
Comment:			

Reagent Analyte Information

Analyte	Source ID	Source Exp. Date	Source Conc.	Source Units	Final Conc.	Final Conc. Units
13C3 HFPO-DA	13C3 HFPO-DA_00004	08/28/2018	50.00000	ug/mL	0.50000	ug/mL
13C3 HFPO-DA (IS)	13C3 HFPO-DA_00004	08/28/2018	50.00000	ug/mL	0.50000	ug/mL

Source Recipients

Reagent	Description	Type	Expiration	Vendor	Vendor Lot #	Vendor Cat Lot #	Volume Used	Volume Units
13C3 HFPO-DA_00004	13C3 HFPO-DA I.S. for HFPO	ASTD	08/28/18	Wellington Laboratories	M3HFPOADA0616M3HFPO-DA	1.00000	mL	

Ok PW
8/29/17

ataset: Untitled

st Altered: Tuesday, August 29, 2017 10:47:21 Mountain Daylight Time

nted: Tuesday, August 29, 2017 10:47:53 Mountain Daylight Time

ethod: C:\MassLynx\8321.PRO\MethDB\hfpo.mdb 23 Aug 2017 10:19:52

libration: C:\MassLynx\8321.PRO\CurveDB\hfpo17d24.cdb 24 Apr 2017 13:20:17

sample Name: hfpo717H23083

PO IS 00004 MRM of 2 channels,ES-
328.8 > 284.8


2.895e+005

HFPO
0.99

27241.82

2.72e4

0.65

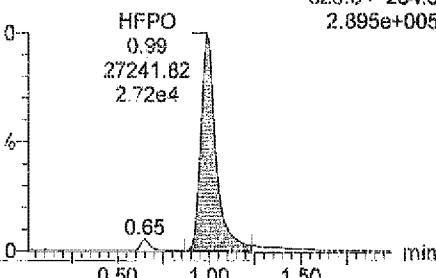


#	Name	Type	Std. Conc.	RT	Area	S/Area	Response	Primar	ppb	%Dev
1	hfpo717H23083		10.000	0.99	27241.822	28795.438	0.946	bd	10.0	-0.4

Dataset: Untitled

Last Altered: Tuesday, August 29, 2017 10:47:21 Mountain Daylight Time
Entered: Tuesday, August 29, 2017 10:47:53 Mountain Daylight Time

Sample Name: hfp0717H23083

PO IS 00004 MRM of 2 channels,ES-
331.8 > 286.8

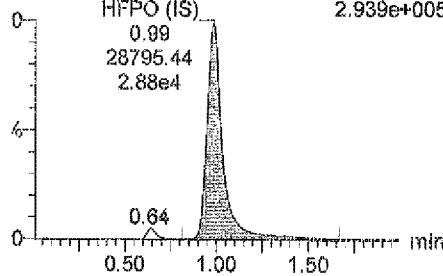

#	Name	Type	Std. Conc.	RT	Area	(S) Area	Response	Primar...	ppb	%Dev
1	hfp0717H23083		1.000	0.99	28795.438		28795.438	bb	1.2	23.6

Reagent

HFPO I.S._00007

Preliminary Report
TestAmerica Denver
Internal Standard Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171212-65681.b\hfpo717L12074.d
 Lims ID: HFPO IS 00007
 Client ID:
 Sample Type: CCV
 Inject. Date: 12-Dec-2017 15:02:32 ALS Bottle#: 25 Worklist Smp#: 74
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: HFPO IS 00007
 Misc. Info.: HFPO17L12
 Operator ID: JBH Instrument ID: LC_LCMS7
 Sublist: chrom-HFPO*sub1
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20171212-65681.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 12-Dec-2017 15:48:38 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d
 Column 1 : Det: F1:MRM
 Process Host: XAWRK024

Averaged ICal Samples:

\\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10026.d
 \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10027.d
 \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10028.d
 \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10029.d
 \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10030.d
 \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10031.d
 \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10032.d
 \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Area Recoveries, Detector: F1:MRM

Compound	Average Standard	Lower Limit	Upper Limit	Sample	% Rec
* 2 13C3 HFPO-DA (IS)	731446	365723	1462892	740105	101.18

RT Recoveries

Compound	Average Standard	Lower Limit	Upper Limit	Sample	DLT(min.)	% Diff
* 2 13C3 HFPO-DA (IS)	0.880	0.380	1.380	1.056	-0.176	19.997

AREA UPPER LIMIT = +100% of internal standard area.

AREA LOWER LIMIT = - 50% of internal standard area.

RT UPPER LIMIT = + 0.500 minutes of internal standard RT.

RT LOWER LIMIT = - 0.500 minutes of internal standard RT.

Reagent

HFPO-DA 00003



WELLINGTON LABORATORIES


CERTIFICATE OF ANALYSIS DOCUMENTATION

PRODUCT CODE:
COMPOUND:

HFPO-DA

2,2,3,3-tetrafluoro-2-(1,1,2,2,3,3-heptafluoropropoxy)-propanoic acid

STRUCTURE: **CAS #:** 13262-13-6

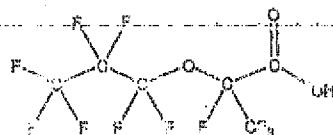

MOLECULAR FORMULA: C₄H₆F₁₀O₃
CONCENTRATION: 50 ± 2.5 µg/ml
CHEMICAL PURITY: >98%
LAST TESTED: (mmddyy) 02/05/2014
EXPIRY DATE: (mmddyy) Stability studies ongoing
RECOMMENDED STORAGE: Store ampoule in a cool, dark place

MOLECULAR WEIGHT: 330.05
SOLVENT(S): Methanol

DOCUMENTATION DATA ATTACHED:

Figure 1: LC/MS Data (TIC and Mass Spectrum)
Figure 2: LC/MS/MS Data (Selected MRM Transitions)

21-0-D25 PB
21-LPL
MDL

ADDITIONAL INFORMATION:

- See page 2 for further details.

FOR LABORATORY USE ONLY: NOT FOR HUMAN OR DRUG USE

Certified By:

B.G. Chittim

Date: 02/13/2014
(mmddyy)

Wellington Laboratories Inc., 346 Southgate Dr. Guelph ON N1G 3M5 CANADA
519-822-2436 • Fax: 519-822-2849 • info@well-labs.com

INTENDED USE:

The products prepared by Wellington Laboratories Inc. are for laboratory use only. They are designed to be used as reference standards for the identification and/or quantification of specific chemical compound(s).

HAZARDS:

This product should only be used by qualified personnel familiar with its potential hazards and trained in the handling of hazardous chemicals. Due care should be exercised to prevent unnecessary human contact or ingestion. All procedures should be carried out in a well-functioning fume hood and suitable gloves, eye protection and clothing should be worn at all times. Waste should be disposed of according to national and regional regulations. Material Safety Data Sheets (MSDSs) are available upon request.

SYNTHESIS / CHARACTERIZATION:

Where possible, all of our products are synthesized using single-product, unambiguous routes. They are then characterized, and their structures and purity confirmed, using a combination of the most relevant techniques, such as NMR, GC/MS, LC/MS/MS, x-ray crystallography and melting point. Isotopic purities of mass-labelled compounds are also confirmed using HRGC/HRMS and/or LC/MS/MS.

HOMOGENEITY:

Prior to solution preparation, crystalline material is tested for homogeneity using a variety of techniques (as stated above) and its solubility in a given solvent is taken into consideration. Duplicate solutions of a new product are prepared from the same crystalline lot and, after the addition of an appropriate internal standard, they are compared by GC/MS and/or LC/MS/MS. The relative response factors of the analyte of interest in each solution are required to be $\pm 5\%$ RSD. New solution lots of existing products are compared to older lots in the same manner, which further confirms the homogeneity of the crystalline material as well as the stability and homogeneity of the solutions in the storage containers.

UNCERTAINTY:

The maximum combined relative standard uncertainty of our reference standard solutions is calculated using the following equation:

The combined relative standard uncertainty, $U_r(V)$, of a value V and the uncertainty of the independent parameters

x_1, x_2, \dots, x_n on which it depends is:

$$U_r(V(x_1, x_2, \dots, x_n)) = \sqrt{\sum_{i=1}^n U(x_i)^2}$$

where U is expressed as a relative standard uncertainty of the individual parameter.

The individual uncertainties taken into account include those associated with weights (calibration of the balance) and volumes (calibration of the volumetric glassware). An expanded maximum combined percent relative uncertainty of $\pm 5\%$ (calculated with a coverage factor of 2 and a level of confidence of 95%) is stated on the Certificate of Analysis for all our products.

TRACEABILITY:

All reference standard solutions are traceable to specific crystalline lots. The microbalances used for solution preparation are regularly tested by an external, ISO/IEC 17020:2006 accredited calibration company, in addition, their calibration is verified prior to each weighing using NIST and/or NIST traceable external weights. All volumetric glassware used is of Class A tolerance and has been tested according to the appropriate ASTM procedures, which are ultimately traceable to NIST. For certain products, traceability to intermediate interlaboratory studies has also been established.

EXPIRY DATE / PERIOD OF VALIDITY:


Ongoing stability studies of this product have demonstrated stability in its composition and concentration for the period of time specified by the expiry date in the unopened ampoule. Monitoring for any degradation or change in concentration of the listed analyte(s) is performed on a routine basis.

LIMITED WARRANTY:

At the time of shipment, all products are warranted to be free of defects in material and workmanship and to conform to the stated technical and purity specifications.



QUALITY MANAGEMENT:

This product was produced using a Quality Management System registered to ISO 9001:2008 by SAI Global, ISO/IEC 17025:2005 by the Canadian Association for Laboratory Accreditation Inc. (CALA, A1226), and ISO GUIDE 34:2009 by ACCLASS (certificate number AF-050).

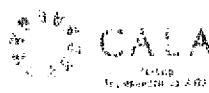

For additional information or assistance concerning this or any other products from Wellington Laboratories Inc., please visit our website at www.well-labs.com or contact us directly at iso@well-labs.com

Figure 1: HFFPO-DA; LC/MS Data (TIC and Mass spectrum)

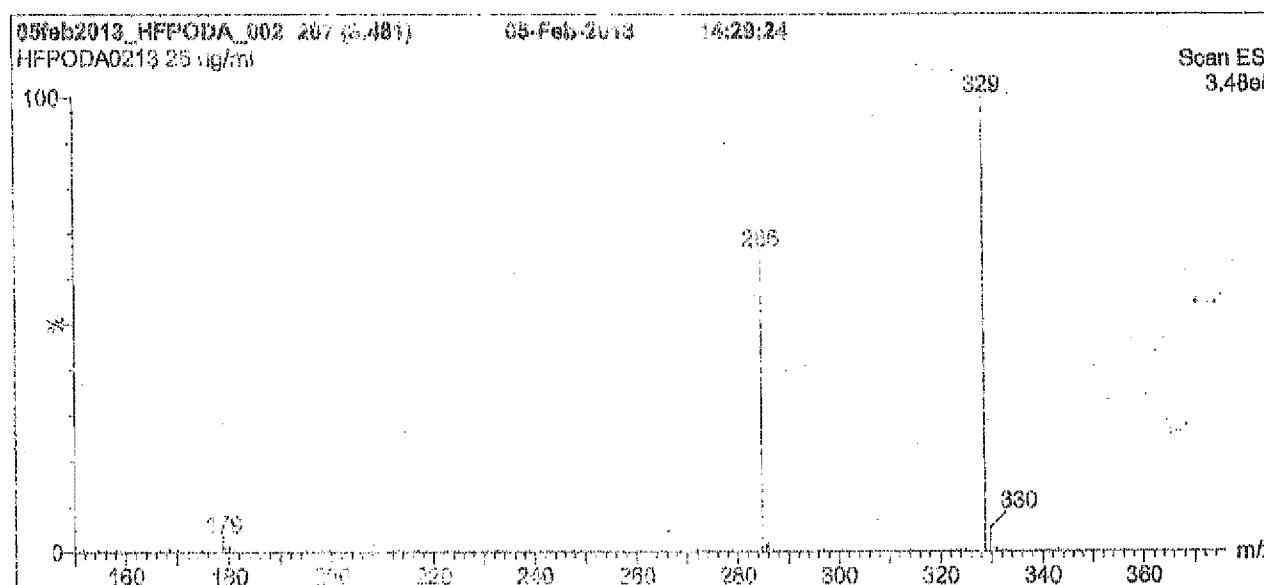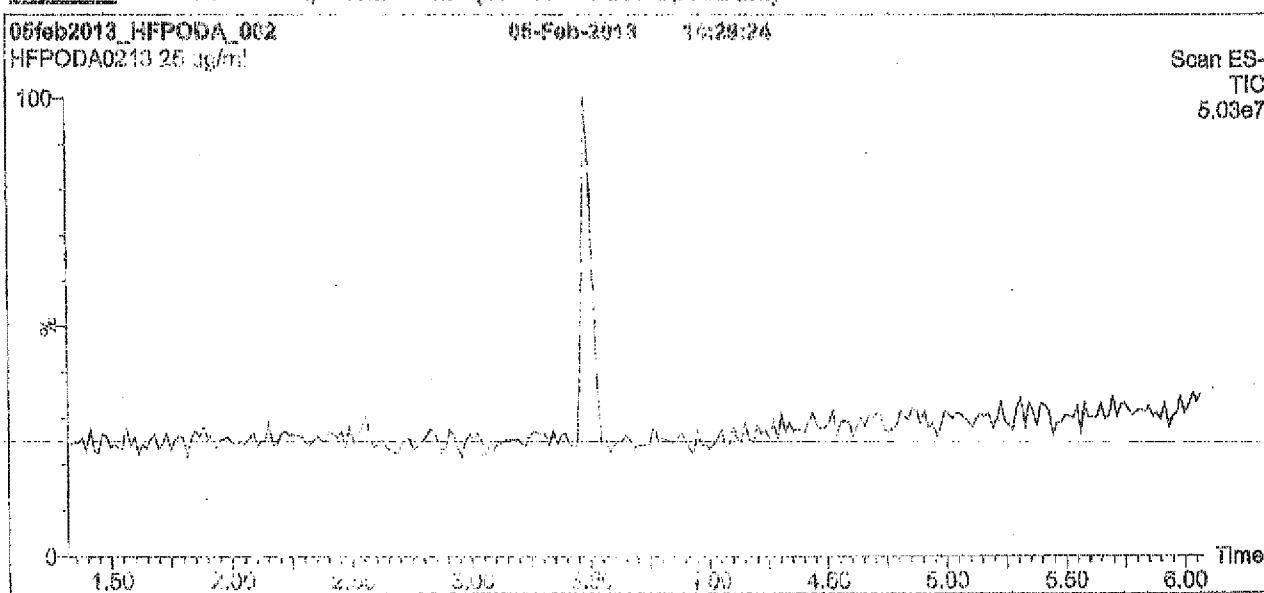

Conditions for Figure 1:

LC: Waters Acuity UHPLC Performance LC
MS: Micromass Quattro micro API MS

Chromatographic Conditions:

Column: Kinetex PEI
2.6 μm , 4.8 x 100 mm

Mobile phase: Gradient
Start: 40% (0.020 Molar LiAc/H₂O / 60% H₂O
Gradient: 10 mM Li₂CO₃ buffer)
Ramp to 60% organic over 6 min and hold for 1 min
before returning to initial conditions in 0.5 min.
Time: 11 min


Flow: 0.01 $\mu\text{l}/\text{min}$

MS Parameters:

Experiment: Full Scan (150 - 850 amu)

Source: Electrospray (negative)
Capillary Voltage (kV) = 3.00
Cone Voltage (V) = 9.00
Cone Gas Flow (l/hr) = 50
Desolvation Gas Flow (l/hr) = 750

Figure 2: HFPO-DA LC/MS/MS Data (Selected MRM Transitions)

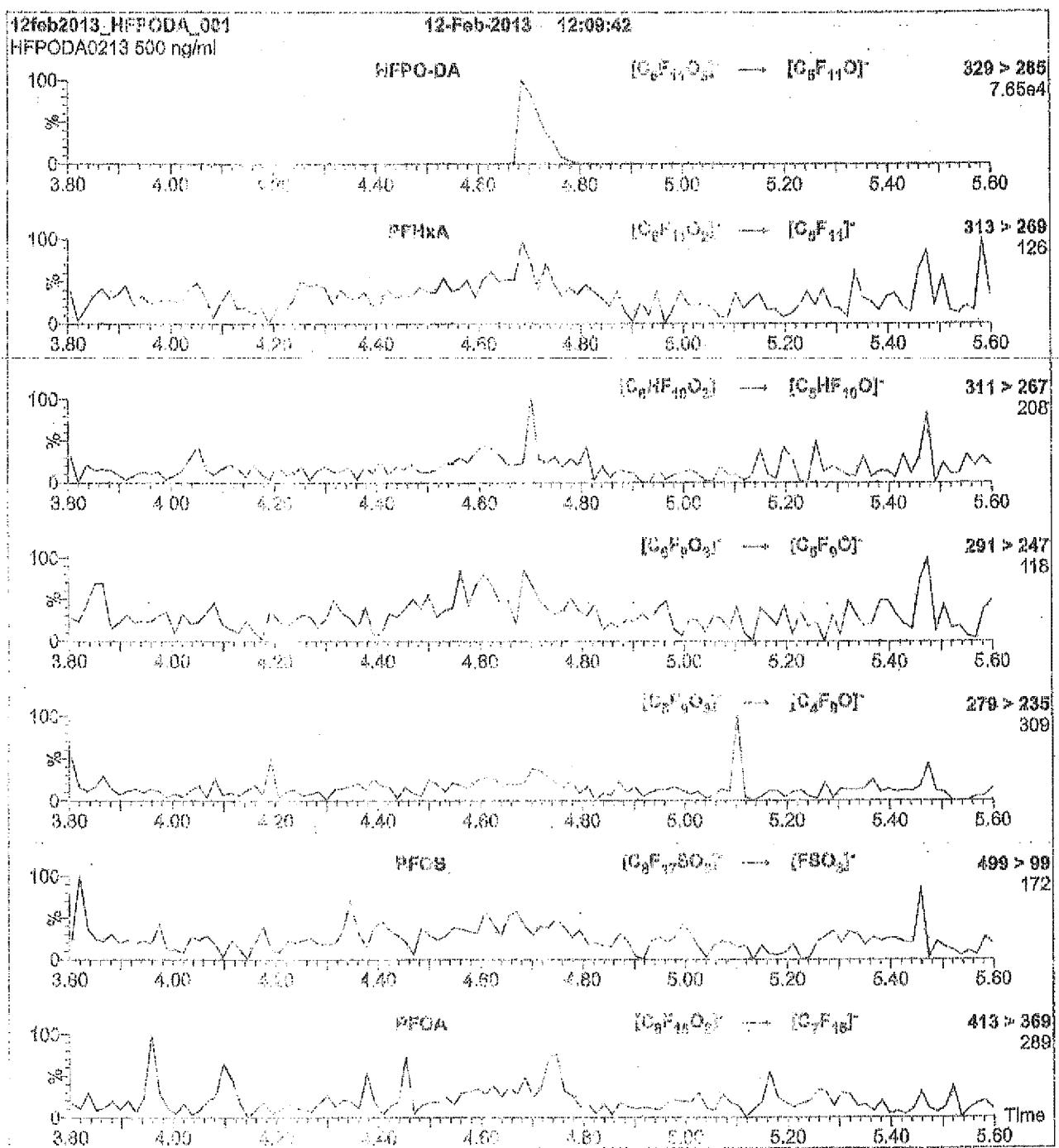

Conditions for Figure 2:

Injection: Direct loop injection
10 μ l (500 ng/ml HFPO-DA)

ESI Parameters:

Collision Gas (mbar) = 3.87e-3

Mobile phase: Isocratic 80% (60:20 MeOH:ACN) / 20% H_2O
(both with 10 mM NH_4OAc buffer)

Collision Energy (eV) = 5

Flow: 300 μ l/min

Reagent

HFPO-DA 00004



WELLINGTON LABORATORIES

CERTIFICATE OF ANALYSIS DOCUMENTATION


PRODUCT CODE: HFPO-DA

LOT NUMBER: HFPODA0717

COMPOUND: 2,3,3,3-Tetrafluoro-2-(1,1,2,2,3,3,3-heptafluoropropoxy)-propanoic acid

STRUCTURE:

CAS #: 13252-13-6

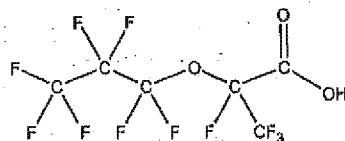

MOLECULAR FORMULA: C₆HF₁₁O₃

MOLECULAR WEIGHT: 330.05

CONCENTRATION: 50 ± 2.5 µg/ml

SOLVENT(S): Methanol

CHEMICAL PURITY: >98%

LAST TESTED: (mm/dd/yyyy) 07/13/2017

EXPIRY DATE: (mm/dd/yyyy) 07/13/2020

RECOMMENDED STORAGE: Store ampoule in a cool, dark place

DOCUMENTATION/ DATA ATTACHED:

Figure 1: LC/MS Data (TIC and Mass Spectrum)

Figure 2: LC/MS/MS Data (Selected MRM Transitions)

ADDITIONAL INFORMATION:

- See page 2 for further details.
- Product is commercially known as GenX.

FOR LABORATORY USE ONLY: NOT FOR HUMAN OR DRUG USE

Certified By:

B.G. Chittim, General Manager

Date: 07/14/2017

(mm/dd/yyyy)

Wellington Laboratories Inc., 345 Southgate Dr. Guelph ON N1G 3M5 CANADA
519-822-2436 • Fax: 519-822-2849 • Info@well-labs.com

INTENDED USE:

The products prepared by Wellington Laboratories Inc. are for laboratory use only. This certified reference material (CRM) was designed to be used as a standard for the identification and/or quantification of the specific chemical compound it contains.

HAZARDS:

This product should only be used by qualified personnel familiar with its potential hazards and trained in the handling of hazardous chemicals. Due care should be exercised to prevent unnecessary human contact or ingestion. All procedures should be carried out in a well-functioning fume hood and suitable gloves, eye protection, and clothing should be worn at all times. Waste should be disposed of according to national and regional regulations. Safety Data Sheets (SDSs) are available upon request.

SYNTHESIS / CHARACTERIZATION:

Where possible, all of our products are synthesized using single-product unambiguous routes. They are then characterized, and their structures and purities confirmed, using a combination of the most relevant techniques, such as NMR, GC/MS, LC/MS/MS, SFC/UV/MS/MS, x-ray crystallography, and melting point. Isotopic purities of mass-labelled compounds are also confirmed using HRGC/HRMS and/or LC/MS/MS.

HOMOGENEITY:

Prior to solution preparation, crystalline material is tested for homogeneity using a variety of techniques (as stated above) and its solubility in a given diluent is taken into consideration. Duplicate solutions of a new product are prepared from the same crystalline lot and, after the addition of an appropriate internal standard, they are compared by GC/MS, LC/MS/MS and/or SFC/UV/MS/MS. The relative response factors of the analyte of interest in each solution are required to be <5% RSD. New solution lots of existing products are compared to older lots in the same manner, which further confirms the homogeneity of the crystalline material as well as the stability and homogeneity of the solutions in the storage containers. In order to maintain the integrity of the assigned value(s), and associated uncertainty, the dilution or injection of a subsample of this product should be performed using calibrated measuring equipment.

UNCERTAINTY:

The maximum combined relative standard uncertainty of our reference standard solutions is calculated using the following equation:

The combined relative standard uncertainty, $u_c(y)$, of a value y and the uncertainty of the independent parameters

x_1, x_2, \dots, x_n on which it depends is:

$$u_c(y(x_1, x_2, \dots, x_n)) = \sqrt{\sum_{i=1}^n u(y, x_i)^2}$$

where x is expressed as a relative standard uncertainty of the individual parameter.

The individual uncertainties taken into account include those associated with weights (calibration of the balance) and volumes (calibration of the volumetric glassware). An expanded maximum combined percent relative uncertainty of ±5% (calculated with a coverage factor of 2 and a level of confidence of 95%) is stated on the Certificate of Analysis for all of our products.

TRACEABILITY:

All reference standard solutions are traceable to specific crystalline lots. The microbalances used for solution preparation are regularly tested by an external ISO/IEC 17025 accredited calibration company. In addition, their calibration is verified prior to each weighing using calibrated NIST and/or NRC traceable external weights. All volumetric glassware used is calibrated, of Class A tolerance, and has been tested according to the appropriate ASTM procedures, which are ultimately traceable to NIST. For certain products, traceability to International Interlaboratory studies has also been established.

EXPIRY DATE / PERIOD OF VALIDITY:

Ongoing stability studies of this product have demonstrated stability in its composition and concentration, until the specified expiry date, in the unopened ampoule. Monitoring for any degradation or change in concentration of the listed analyte(s) is performed on a routine basis.

LIMITED WARRANTY:

At the time of shipment, all products are warranted to be free of defects in material and workmanship and to conform to the stated technical and purity specifications.



QUALITY MANAGEMENT:

This product was produced using a Quality Management System registered to the latest versions of ISO 9001 by SAI Global, ISO/IEC 17025 by the Canadian Association for Laboratory Accreditation Inc. (CALA; A 1226), and ISO GUIDE 34 by ANSI-ASQ National Accreditation Board (ANAB; AR-1523).



For additional information or assistance concerning this or any other products from Wellington Laboratories Inc., please visit our website at www.well-labs.com or contact us directly at info@well-labs.com

Figure 1: HFPO-DA; LC/MS Data (TIC and Mass Spectrum)

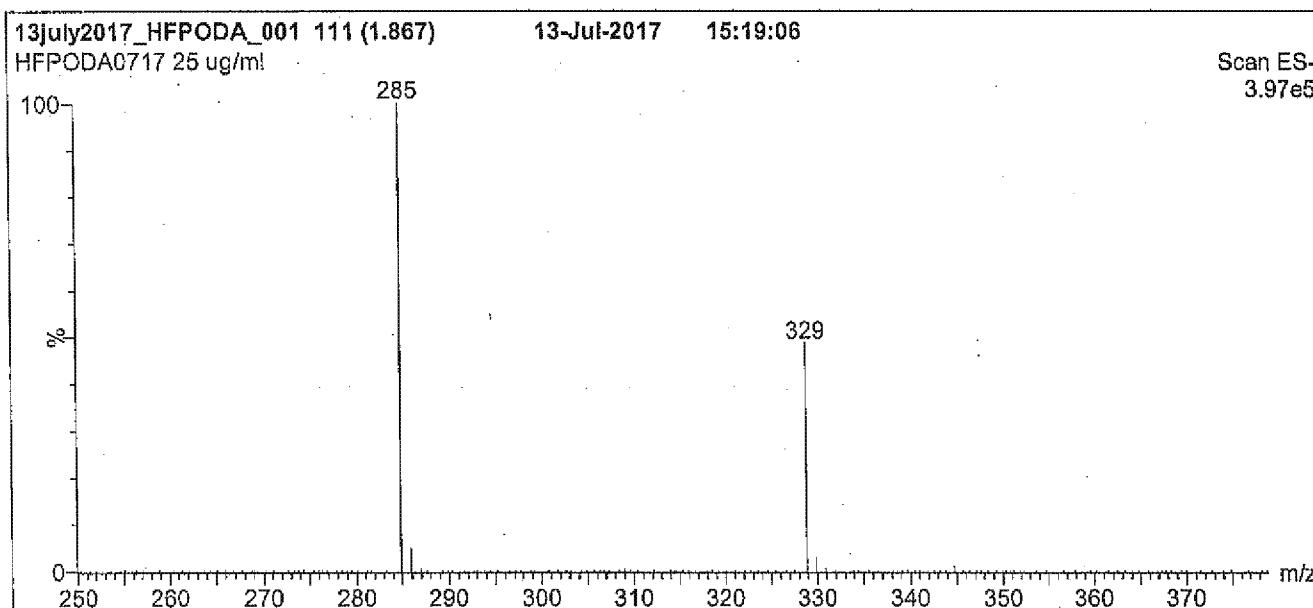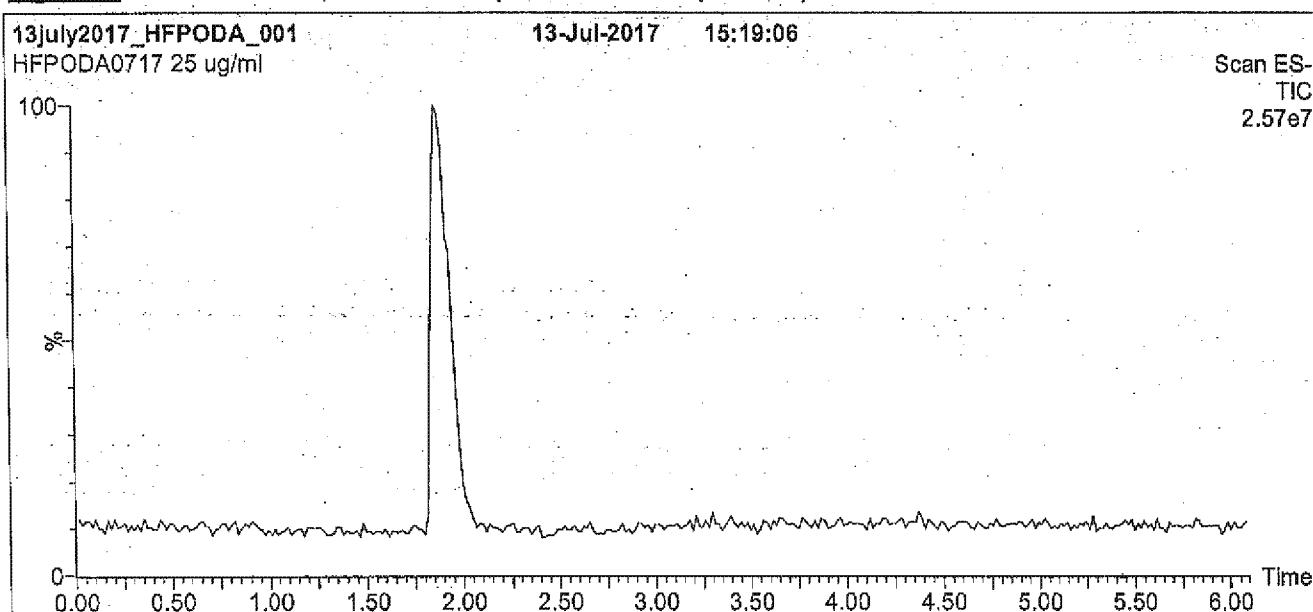

Conditions for Figure 1:

LC: Waters Acquity Ultra Performance LC
MS: Micromass Quattro *micro* API MS

Chromatographic Conditions

Column: Acquity UPLC BEH Shield RP₁₈
1.7 μ m, 2.1 x 100 mm

Mobile phase: Gradient
Start: 55% MeOH / 45% H₂O with 10 mM NH₄OAc buffer
Ramp to 90% organic over 7.5 min and hold for 1.5 min
before returning to initial conditions in 0.5 min.

Time: 10 min


Flow: 300 μ l/min

MS Parameters

Experiment: Full Scan (250 - 850 amu)

Source: Electrospray (negative)
Capillary Voltage (kV) = 3.00
Cone Voltage (V) = 10.00
Cone Gas Flow (l/hr) = 100
Desolvation Gas Flow (l/hr) = 700

Figure 2: HFPO-DA; LC/MS/MS Data (Selected MRM Transitions)

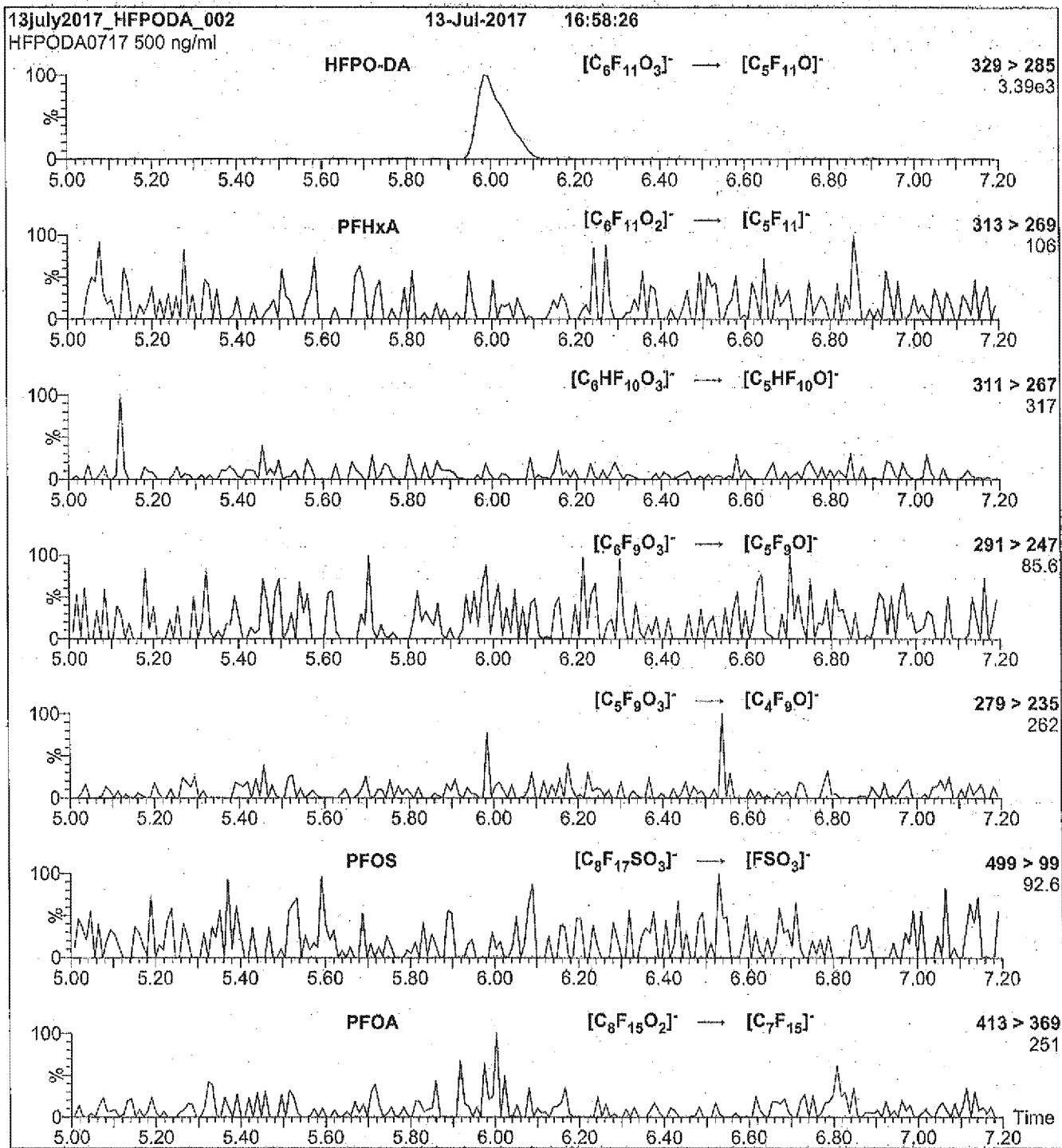

Conditions for Figure 2:

Injection: Direct loop injection
10 μ l (500 ng/ml HFPO-DA)

MS Parameters

Collision Gas (mbar) = 3.20e-3
Collision Energy (eV) = 5

Mobile phase: Isocratic 80% MeOH / 20% H₂O with 10 mM NH₄OAc buffer

Flow: 300 μ l/min

8321A_HFPO_Du

HFPO-DA

FORM II
LCMS SURROGATE RECOVERY

Lab Name: TestAmerica Denver

Job No.: 320-35355-1

SDG No.: _____

Matrix: Water

Level: Low

GC Column (1): Synergi Hyd ID: _____

Client Sample ID	Lab Sample ID	HFPEDA #
FAY-D-7154BUTLE-W1 -1-012418	320-35355-1	103
FAY-D-384JERML-W1- 1-012418	320-35355-2	107
FAY-D-2528INDIA-W1- -2-012418	320-35355-3	104
FAY-D-661PINEA-W1- 1-012418	320-35355-4	97
FAY-D-701PINEA-W1- 1-012418	320-35355-5	102
FAY-D-701PINEA-W1- 2-012418	320-35355-6	125
FAY-D-663PINEA-W1- 1-012418	320-35355-7	123
FAY-D-7166BUTLE-W2 -1-012418	320-35355-8	100
FAY-D-7166BUTLE-W1 -1-012418	320-35355-9	90
FAY-D-7166BUTLE-W1 -1-012418D	320-35355-10	103
	MB 280-403139/1-A	111
	LCS 280-403139/2-A	107
	LCSD 280-403139/3-A	105
	LLCS 280-403139/4-A	128
FAY-D-7166BUTLE-W1 -1-012418 MS	320-35355-9 MS	123
FAY-D-7166BUTLE-W1 -1-012418 DU	320-35355-9 DU	100
	DLCK 280-390728/12	102

HFPEDA = 13C3 HFPO-DA

QC LIMITS
50-200

Column to be used to flag recovery values

FORM II 8321A

FORM III
LCMS LAB CONTROL SAMPLE RECOVERY

Lab Name: TestAmerica Denver Job No.: 320-35355-1

SDG No.: _____

Matrix: Water Level: Low Lab File ID: hfpo718A30005.d

Lab ID: LCS 280-403139/2-A Client ID: _____

COMPOUND	SPIKE ADDED (ug/L)	LCS CONCENTRATION (ug/L)	LCS % REC	QC LIMITS REC	#
HFPO-DA	0.200	0.225	113	70-130	

Column to be used to flag recovery and RPD values

FORM III 8321A

FORM III
LCMS LAB CONTROL SAMPLE DUPLICATE RECOVERY

Lab Name: TestAmerica Denver

Job No.: 320-35355-1

SDG No.: _____

Matrix: Water Level: Low Lab File ID: hfpo718A30006.d

Lab ID: LCSD 280-403139/3-A Client ID: _____

COMPOUND	SPIKE ADDED (ug/L)	LCSD CONCENTRATION (ug/L)	LCSD REC	%	QC LIMITS		#
					RPD	RPD	
HFPO-DA	0.200	0.207	103	8	20	70-130	

Column to be used to flag recovery and RPD values

FORM III 8321A

FORM III
LCMS LOW LEVEL CONTROL SAMPLE RECOVERY

Lab Name: TestAmerica Denver

Job No.: 320-35355-1

SDG No.: _____

Matrix: Water Level: Low Lab File ID: hfpo718A30007.d

Lab ID: LLCS 280-403139/4-A Client ID: _____

COMPOUND	SPIKE ADDED (ug/L)	LLCS CONCENTRATION (ug/L)	LLCS REC	QC LIMITS REC	#
HFPO-DA	0.0200	0.0147	74	70-130	

Column to be used to flag recovery and RPD values

FORM III 8321A

FORM III
LCMS MATRIX SPIKE RECOVERY

Lab Name: TestAmerica Denver Job No.: 320-35355-1

SDG No.: _____

Matrix: Water Level: Low Lab File ID: hfpo718A30019.d

Lab ID: 320-35355-9 MS Client ID: FAY-D-7166BUTLE-W1-1-012418 MS

COMPOUND	SPIKE ADDED (ug/L)	SAMPLE CONCENTRATION (ug/L)	MS CONCENTRATION (ug/L)	MS % REC	QC LIMITS REC	#
HFPO-DA	0.203	0.017	0.202	91	70-130	

Column to be used to flag recovery and RPD values

FORM III 8321A

FORM III
LCMS DETECTION LIMIT CHECK STANDARD RECOVERY

Lab Name: TestAmerica Denver Job No.: 320-35355-1

SDG No.: _____

Matrix: Water Level: Low Lab File ID: hfpo717J10035.d

Lab ID: DLCK 280-390728/12 Client ID: _____

COMPOUND	SPIKE ADDED (ug/L)	DLCK CONCENTRATION (ug/L)	DLCK % REC	QC LIMITS REC	#
HFPO-DA	0.250	<0.50	78	70-130	

Column to be used to flag recovery and RPD values

FORM III 8321A

FORM IV
LCMS METHOD BLANK SUMMARY

Lab Name: TestAmerica Denver Job No.: 320-35355-1
SDG No.:
Lab File ID: hfpo718A30004.d Lab Sample ID: MB 280-403139/1-A
Matrix: Water Date Extracted: 01/29/2018 21:37
Instrument ID: LC_LCMS7 Date Analyzed: 01/30/2018 09:15
Level: (Low/Med) Low

THIS METHOD BLANK APPLIES TO THE FOLLOWING SAMPLES:

CLIENT SAMPLE ID	LAB SAMPLE ID	LAB FILE ID	DATE ANALYZED
	LCS 280-403139/2-A	hfpo718A300 05.d	01/30/2018 09:19
	LCSD 280-403139/3-A	hfpo718A300 06.d	01/30/2018 09:22
	LLCS 280-403139/4-A	hfpo718A300 07.d	01/30/2018 09:25
FAY-D-7154BUTLE-W1-1-012418	320-35355-1	hfpo718A300 08.d	01/30/2018 09:28
FAY-D-384JERML-W1-1-012418	320-35355-2	hfpo718A300 09.d	01/30/2018 09:32
FAY-D-2528INDIA-W1-2-012418	320-35355-3	hfpo718A300 10.d	01/30/2018 09:35
FAY-D-661PINEA-W1-1-012418	320-35355-4	hfpo718A300 11.d	01/30/2018 09:38
FAY-D-701PINEA-W1-1-012418	320-35355-5	hfpo718A300 13.d	01/30/2018 09:45
FAY-D-701PINEA-W1-2-012418	320-35355-6	hfpo718A300 14.d	01/30/2018 09:48
FAY-D-663PINEA-W1-1-012418	320-35355-7	hfpo718A300 15.d	01/30/2018 09:51
FAY-D-7166BUTLE-W2-1-012418	320-35355-8	hfpo718A300 16.d	01/30/2018 09:54
FAY-D-7166BUTLE-W1-1-012418	320-35355-9	hfpo718A300 17.d	01/30/2018 09:58
FAY-D-7166BUTLE-W1-1-012418 DU	320-35355-9 DU	hfpo718A300 18.d	01/30/2018 10:01
FAY-D-7166BUTLE-W1-1-012418 MS	320-35355-9 MS	hfpo718A300 19.d	01/30/2018 10:04
FAY-D-7166BUTLE-W1-1-012418D	320-35355-10	hfpo718A300 20.d	01/30/2018 10:07

FORM I
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: TestAmerica Denver Job No.: 320-35355-1

SDG No.: _____

Client Sample ID: FAY-D-7154BUTLE-W1-1-0124 Lab Sample ID: 320-35355-1
18 _____

Matrix: Water Lab File ID: hfpo718A30008.d

Analysis Method: 8321A Date Collected: 01/24/2018 15:19

Extraction Method: 3535 Date Extracted: 01/29/2018 21:37

Sample wt/vol: 251.6 (mL) Date Analyzed: 01/30/2018 09:28

Con. Extract Vol.: 5 (mL) Dilution Factor: 1

Injection Volume: 20 (uL) GC Column: Synergi Hydro ID: _____

% Moisture: GPC Cleanup: (Y/N) N

Analysis Batch No.: 403246 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	0.047		0.010	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	103		50-200

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30008.d
 Lims ID: 320-35355-B-1-A
 Client ID: FAY-D-7154BUTLE-W1-1-012418
 Sample Type: Client
 Inject. Date: 30-Jan-2018 09:28:48 ALS Bottle#: 15 Worklist Smp#: 8
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: 320-35355-B-1-A
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 30-Jan-2018 14:49:56 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK010

First Level Reviewer: meyera Date: 30-Jan-2018 14:47:46

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA
 331.8 > 286.8 0.920 0.961 -0.041 1.000 750111 10.3 1334
 * 2 13C3 HFPO-DA (IS)
 331.8 > 286.8 0.920 0.961 -0.041 750111 10.0 1334
 1 Perfluoro(2-propoxypropanoic) acid
 328.8 > 284.8 0.934 0.988 -0.054 1.000 196846 2.38 47.9


TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20180130-66815.b\\hfpo718A30008.d
Injection Date: 30-Jan-2018 09:28:48 Instrument ID: LC_LCMS7
Lims ID: 320-35355-B-1-A Lab Sample ID: 280-35355-1
Client ID: FAY-D-7154BUTLE-W1-1-012418
Operator ID: JBH ALS Bottle#: 15 Worklist Smp#: 8
Injection Vol: 20.0 ul Dil. Factor: 1.0000
Method: HFPO Limit Group: LC - 8321A_HFPO_Du

\$ 3 13C3 HFPO-DA

* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

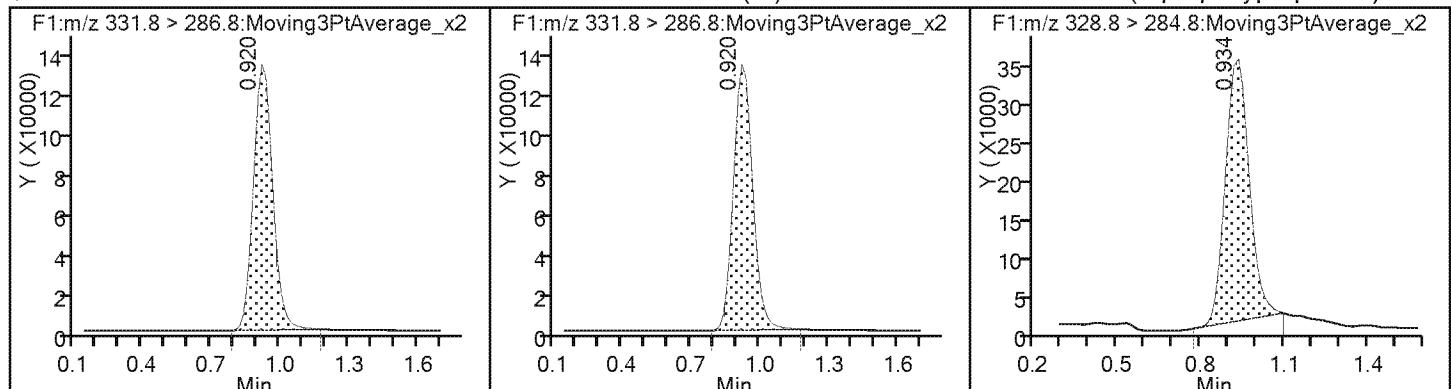

TestAmerica Denver
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30008.d
 Lims ID: 320-35355-B-1-A
 Client ID: FAY-D-7154BUTLE-W1-1-012418
 Sample Type: Client
 Inject. Date: 30-Jan-2018 09:28:48 ALS Bottle#: 15 Worklist Smp#: 8
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: 320-35355-B-1-A
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 30-Jan-2018 14:49:56 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK010

First Level Reviewer: meyera Date: 30-Jan-2018 14:47:46

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	10.3	102.55

FORM I
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: TestAmerica Denver Job No.: 320-35355-1

SDG No.: _____

Client Sample ID: FAY-D-384JERML-W1-1-01241 Lab Sample ID: 320-35355-2
8

Matrix: Water Lab File ID: hfpo718A30009.d

Analysis Method: 8321A Date Collected: 01/24/2018 13:35

Extraction Method: 3535 Date Extracted: 01/29/2018 21:37

Sample wt/vol: 253 (mL) Date Analyzed: 01/30/2018 09:32

Con. Extract Vol.: 5 (mL) Dilution Factor: 1

Injection Volume: 20 (uL) GC Column: Synergi Hydro ID: _____

% Moisture: GPC Cleanup: (Y/N) N

Analysis Batch No.: 403246 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	<0.010		0.010	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	107		50-200

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30009.d
 Lims ID: 320-35355-B-2-A
 Client ID: FAY-D-384JERML-W1-1-012418
 Sample Type: Client
 Inject. Date: 30-Jan-2018 09:32:03 ALS Bottle#: 16 Worklist Smp#: 9
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: 320-35355-B-2-A
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 30-Jan-2018 14:49:56 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK010


First Level Reviewer: meyera Date: 30-Jan-2018 14:47:51

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA
 331.8 > 286.8 0.920 0.961 -0.041 1.000 779773 10.7 1590
 * 2 13C3 HFPO-DA (IS)
 331.8 > 286.8 0.920 0.961 -0.041 1.000 779773 10.0 1590
 1 Perfluoro(2-propoxypropanoic) acid
 328.8 > 284.8 0.934 0.988 -0.054 1.000 40739 0.3003 11.2

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20180130-66815.b\\hfpo718A30009.d
Injection Date: 30-Jan-2018 09:32:03 Instrument ID: LC_LCMS7
Lims ID: 320-35355-B-2-A Lab Sample ID: 280-35355-2
Client ID: FAY-D-384JERML-W1-1-012418
Operator ID: JBH ALS Bottle#: 16 Worklist Smp#: 9
Injection Vol: 20.0 ul Dil. Factor: 1.0000
Method: HFPO Limit Group: LC - 8321A_HFPO_Du

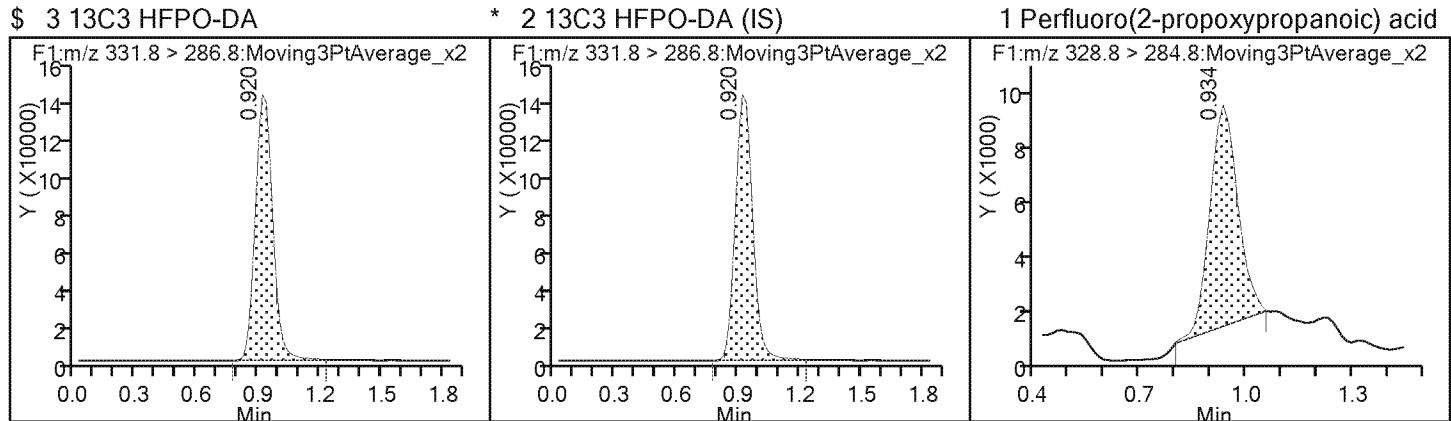

TestAmerica Denver
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30009.d
 Lims ID: 320-35355-B-2-A
 Client ID: FAY-D-384JERML-W1-1-012418
 Sample Type: Client
 Inject. Date: 30-Jan-2018 09:32:03 ALS Bottle#: 16 Worklist Smp#: 9
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: 320-35355-B-2-A
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 30-Jan-2018 14:49:56 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK010

First Level Reviewer: meyera Date: 30-Jan-2018 14:47:51

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	10.7	106.61

FORM I
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: TestAmerica Denver Job No.: 320-35355-1

SDG No.: _____

Client Sample ID: FAY-D-2528INDIA-W1-2-0124 Lab Sample ID: 320-35355-3
18

Matrix: Water Lab File ID: hfpo718A30010.d

Analysis Method: 8321A Date Collected: 01/24/2018 08:56

Extraction Method: 3535 Date Extracted: 01/29/2018 21:37

Sample wt/vol: 251.3 (mL) Date Analyzed: 01/30/2018 09:35

Con. Extract Vol.: 5 (mL) Dilution Factor: 1

Injection Volume: 20 (uL) GC Column: Synergi Hydro ID: _____

% Moisture: GPC Cleanup: (Y/N) N

Analysis Batch No.: 403246 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	<0.010		0.010	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	104		50-200

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30010.d
 Lims ID: 320-35355-B-3-A
 Client ID: FAY-D-2528INDIA-W1-2-012418
 Sample Type: Client
 Inject. Date: 30-Jan-2018 09:35:17 ALS Bottle#: 17 Worklist Smp#: 10
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: 320-35355-B-3-A
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 30-Jan-2018 14:49:56 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK010

First Level Reviewer: meyera Date: 30-Jan-2018 14:47:55

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA
 331.8 > 286.8 0.920 0.961 -0.041 1.000 763513 10.4 1082
 * 2 13C3 HFPO-DA (IS)
 331.8 > 286.8 0.920 0.961 -0.041 1.000 763513 10.0 1082
 1 Perfluoro(2-propoxypropanoic) acid
 328.8 > 284.8 0.934 0.988 -0.054 1.000 30228 0.1753 11.4


TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20180130-66815.b\\hfpo718A30010.d
Injection Date: 30-Jan-2018 09:35:17 Instrument ID: LC_LCMS7
Lims ID: 320-35355-B-3-A Lab Sample ID: 280-35355-3
Client ID: FAY-D-2528INDIA-W1-2-012418
Operator ID: JBH ALS Bottle#: 17 Worklist Smp#: 10
Injection Vol: 20.0 ul Dil. Factor: 1.0000
Method: HFPO Limit Group: LC - 8321A_HFPO_Du

\$ 3 13C3 HFPO-DA

* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

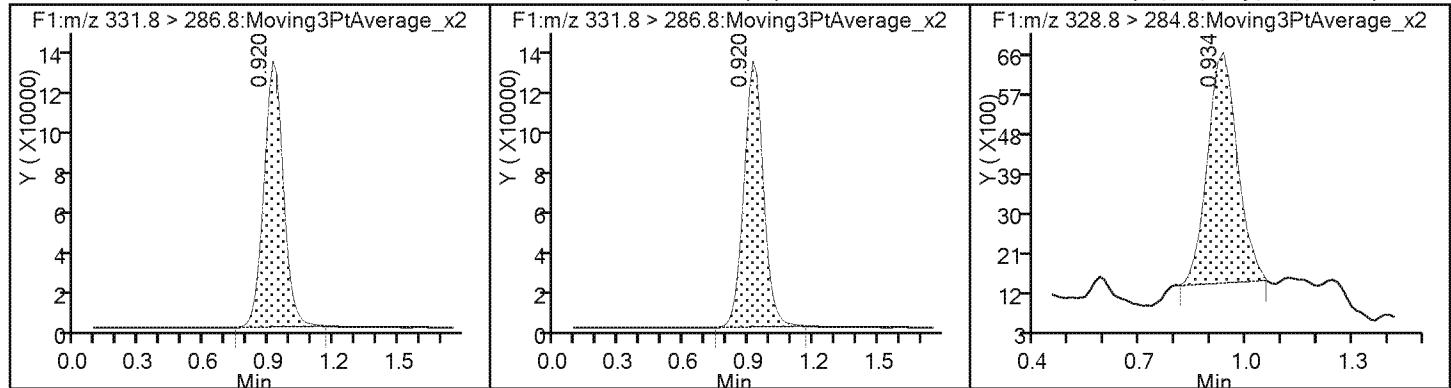

TestAmerica Denver
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30010.d
 Lims ID: 320-35355-B-3-A
 Client ID: FAY-D-2528INDIA-W1-2-012418
 Sample Type: Client
 Inject. Date: 30-Jan-2018 09:35:17 ALS Bottle#: 17 Worklist Smp#: 10
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: 320-35355-B-3-A
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 30-Jan-2018 14:49:56 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK010

First Level Reviewer: meyera Date: 30-Jan-2018 14:47:55

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	10.4	104.38

FORM I
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: TestAmerica Denver Job No.: 320-35355-1

SDG No.:

Client Sample ID: FAY-D-661PINEA-W1-1-01241 Lab Sample ID: 320-35355-4
8

Matrix: Water Lab File ID: hfpo718A30011.d

Analysis Method: 8321A Date Collected: 01/24/2018 09:40

Extraction Method: 3535 Date Extracted: 01/29/2018 21:37

Sample wt/vol: 236.4 (mL) Date Analyzed: 01/30/2018 09:38

Con. Extract Vol.: 5 (mL) Dilution Factor: 1

Injection Volume: 20 (uL) GC Column: Synergi Hydro ID:

% Moisture: GPC Cleanup: (Y/N) N

Analysis Batch No.: 403246 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	0.040		0.010	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	97		50-200

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30011.d
 Lims ID: 320-35355-B-4-A
 Client ID: FAY-D-661PINEA-W1-1-012418
 Sample Type: Client
 Inject. Date: 30-Jan-2018 09:38:32 ALS Bottle#: 18 Worklist Smp#: 11
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: 320-35355-B-4-A
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 30-Jan-2018 14:49:56 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK010


First Level Reviewer: meyera Date: 30-Jan-2018 14:47:58

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA
 331.8 > 286.8 0.920 0.961 -0.041 1.000 707573 9.67 1250
 * 2 13C3 HFPO-DA (IS)
 331.8 > 286.8 0.920 0.961 -0.041 1.000 707573 10.0 1250
 1 Perfluoro(2-propoxypropanoic) acid
 328.8 > 284.8 0.934 0.988 -0.054 1.000 150941 1.89 60.9

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20180130-66815.b\\hfpo718A30011.d
Injection Date: 30-Jan-2018 09:38:32 Instrument ID: LC_LCMS7
Lims ID: 320-35355-B-4-A Lab Sample ID: 280-35355-4
Client ID: FAY-D-661PINEA-W1-1-012418
Operator ID: JBH ALS Bottle#: 18 Worklist Smp#: 11
Injection Vol: 20.0 ul Dil. Factor: 1.0000
Method: HFPO Limit Group: LC - 8321A_HFPO_Du

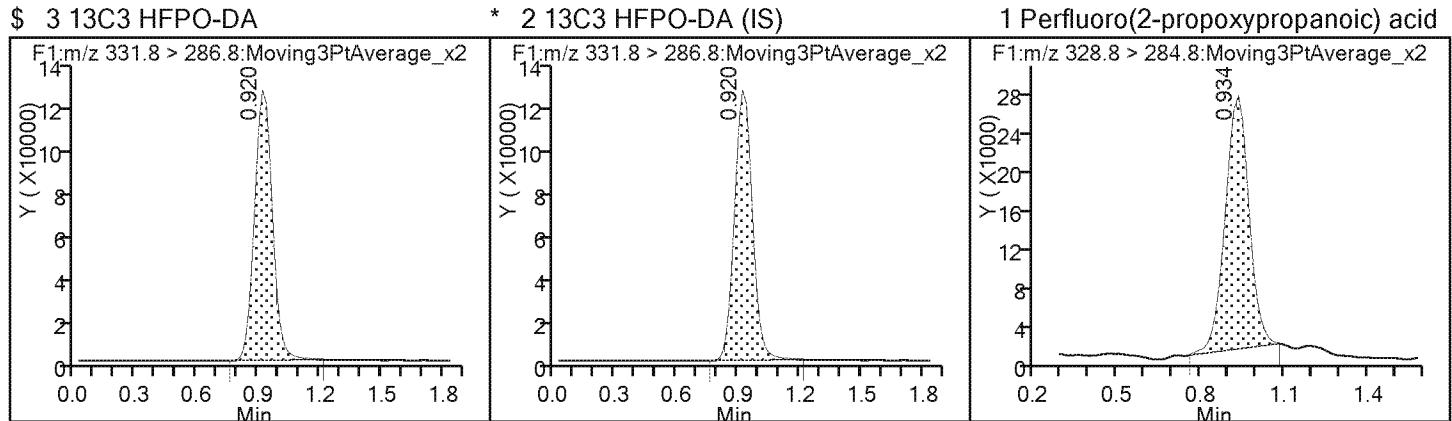

TestAmerica Denver
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30011.d
 Lims ID: 320-35355-B-4-A
 Client ID: FAY-D-661PINEA-W1-1-012418
 Sample Type: Client
 Inject. Date: 30-Jan-2018 09:38:32 ALS Bottle#: 18 Worklist Smp#: 11
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: 320-35355-B-4-A
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 30-Jan-2018 14:49:56 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK010

First Level Reviewer: meyera Date: 30-Jan-2018 14:47:58

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	9.67	96.74

FORM I
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: TestAmerica Denver Job No.: 320-35355-1

SDG No.: _____

Client Sample ID: FAY-D-701PINEA-W1-1-01241 Lab Sample ID: 320-35355-5
8

Matrix: Water Lab File ID: hfpo718A30013.d

Analysis Method: 8321A Date Collected: 01/24/2018 09:58

Extraction Method: 3535 Date Extracted: 01/29/2018 21:37

Sample wt/vol: 240.5 (mL) Date Analyzed: 01/30/2018 09:45

Con. Extract Vol.: 5 (mL) Dilution Factor: 1

Injection Volume: 20 (uL) GC Column: Synergi Hydro ID: _____

% Moisture: GPC Cleanup: (Y/N) N

Analysis Batch No.: 403246 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	<0.010		0.010	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	102		50-200

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30013.d
 Lims ID: 320-35355-B-5-A
 Client ID: FAY-D-701PINEA-W1-1-012418
 Sample Type: Client
 Inject. Date: 30-Jan-2018 09:45:02 ALS Bottle#: 19 Worklist Smp#: 13
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: 320-35355-B-5-A
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 30-Jan-2018 14:50:04 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK010

First Level Reviewer: meyera Date: 30-Jan-2018 14:48:25

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 > 286.8 0.961 0.961 0.0 1.000 745632 10.2 773

* 2 13C3 HFPO-DA (IS)

331.8 > 286.8 0.961 0.961 0.0 1.000 745632 10.0 773

1 Perfluoro(2-propoxypropanoic) acid M

328.8 > 284.8 0.974 0.988 -0.014 1.000 26929 0.1410 5.3 M

QC Flag Legend

Review Flags

M - Manually Integrated

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20180130-66815.b\\hfp0718A30013.d

Injection Date: 30-Jan-2018 09:45:02

Instrument ID: LC_LCMS7

Lims ID: 320-35355-B-5-A

Lab Sample ID: 280-35355-5

Client ID: FAY-D-701PINEA-W1-1-012418

Operator ID: JBH

ALS Bottle#: 19 Worklist Smp#: 13

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A_HFPO_Du

\$ 3 13C3 HFPO-DA

* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid (M)

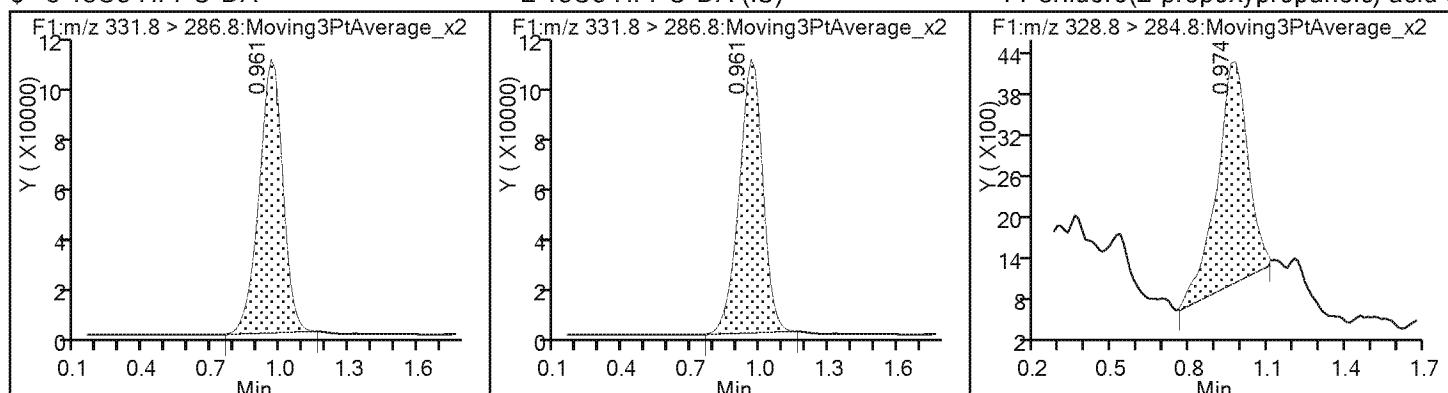

TestAmerica Denver
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30013.d
 Lims ID: 320-35355-B-5-A
 Client ID: FAY-D-701PINEA-W1-1-012418
 Sample Type: Client
 Inject. Date: 30-Jan-2018 09:45:02 ALS Bottle#: 19 Worklist Smp#: 13
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: 320-35355-B-5-A
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 30-Jan-2018 14:50:04 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK010

First Level Reviewer: meyera Date: 30-Jan-2018 14:48:25

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	10.2	101.94


TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20180130-66815.b\\hfpo718A30013.d
 Injection Date: 30-Jan-2018 09:45:02 Instrument ID: LC_LCMS7
 Lims ID: 320-35355-B-5-A Lab Sample ID: 280-35355-5
 Client ID: FAY-D-701PINEA-W1-1-012418
 Operator ID: JBH ALS Bottle#: 19 Worklist Smp#: 13
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Method: HFPO Limit Group: LC - 8321A_HFPO_Du
 Column: Detector F1:MRM

1 Perfluoro(2-propoxypropanoic) acid, CAS: 13252-13-6
 Signal: 1


Not Detected
 Expected RT: 0.99

Processing Integration Results

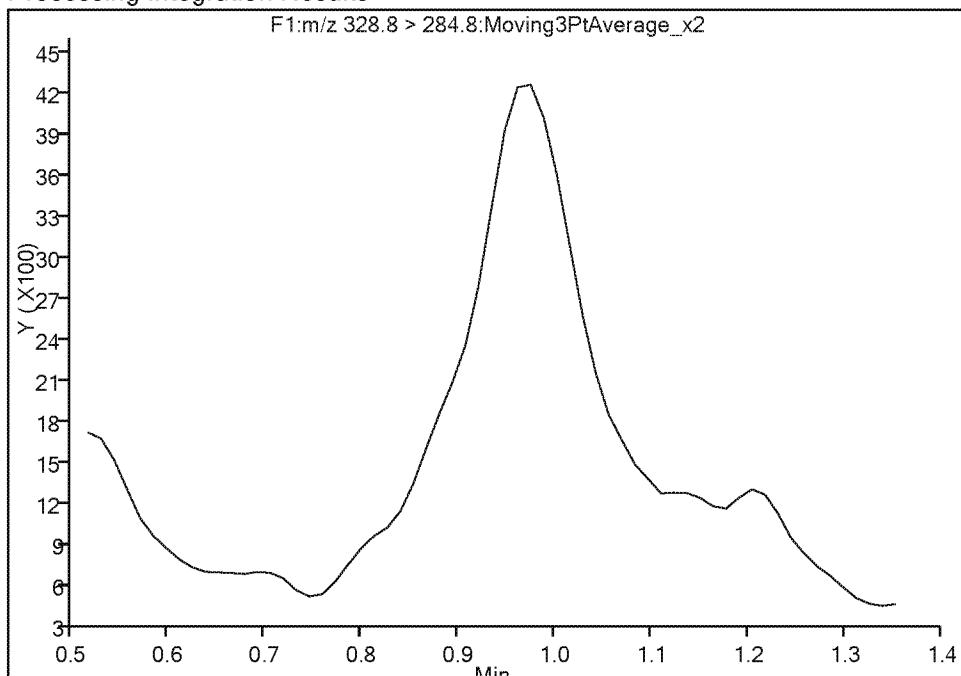

RT: 0.97
 Area: 26929
 Amount: 0.140971
 Amount Units: ug/l

Manual Integration Results

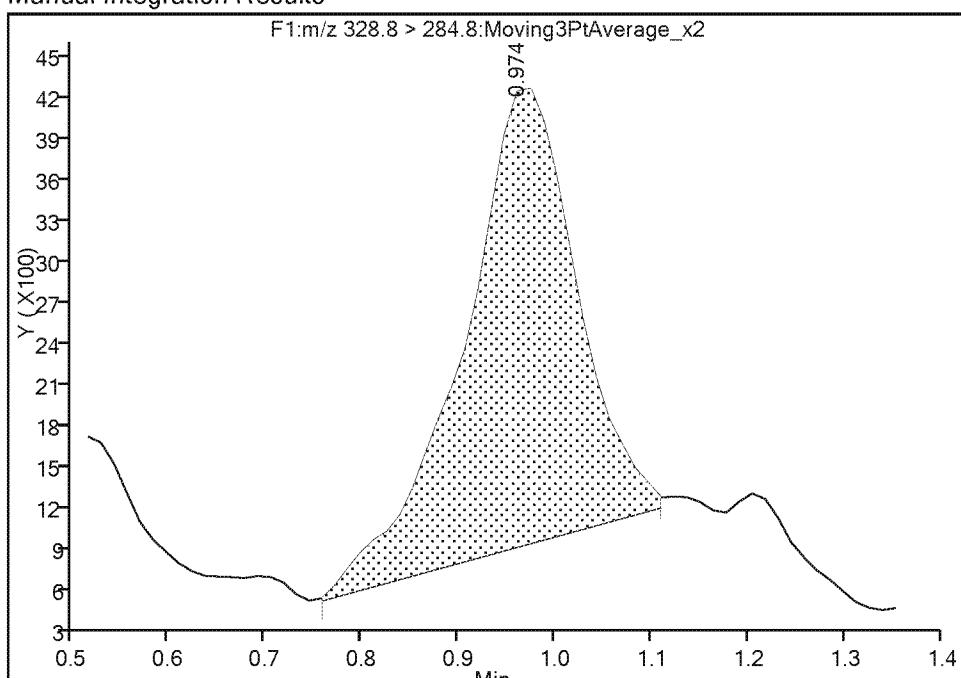

Reviewer: meyera, 30-Jan-2018 14:48:23

Audit Action: Manually Integrated

Audit Reason: Assign Peak

FORM I
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: TestAmerica Denver Job No.: 320-35355-1

SDG No.: _____

Client Sample ID: FAY-D-701PINEA-W1-2-01241 Lab Sample ID: 320-35355-6
8

Matrix: Water Lab File ID: hfpo718A30014.d

Analysis Method: 8321A Date Collected: 01/24/2018 10:05

Extraction Method: 3535 Date Extracted: 01/29/2018 21:37

Sample wt/vol: 252.3 (mL) Date Analyzed: 01/30/2018 09:48

Con. Extract Vol.: 5 (mL) Dilution Factor: 1

Injection Volume: 20 (uL) GC Column: Synergi Hydro ID: _____

% Moisture: GPC Cleanup: (Y/N) N

Analysis Batch No.: 403246 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	<0.010		0.010	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	125		50-200

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30014.d
 Lims ID: 320-35355-B-6-A
 Client ID: FAY-D-701PINEA-W1-2-012418
 Sample Type: Client
 Inject. Date: 30-Jan-2018 09:48:17 ALS Bottle#: 20 Worklist Smp#: 14
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: 320-35355-B-6-A
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 30-Jan-2018 14:50:04 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK010

First Level Reviewer: meyera Date: 30-Jan-2018 14:48:31

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA
 331.8 > 286.8 0.920 0.961 -0.041 1.000 917688 12.5 1638
 * 2 13C3 HFPO-DA (IS)
 331.8 > 286.8 0.920 0.961 -0.041 917688 10.0 1638


TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20180130-66815.b\\hfpo718A30014.d
Injection Date: 30-Jan-2018 09:48:17 Instrument ID: LC_LCMS7
Lims ID: 320-35355-B-6-A Lab Sample ID: 280-35355-6
Client ID: FAY-D-701PINEA-W1-2-012418
Operator ID: JBH ALS Bottle#: 20 Worklist Smp#: 14
Injection Vol: 20.0 ul Dil. Factor: 1.0000
Method: HFPO Limit Group: LC - 8321A_HFPO_Du

\$ 3 13C3 HFPO-DA

* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid (ND)

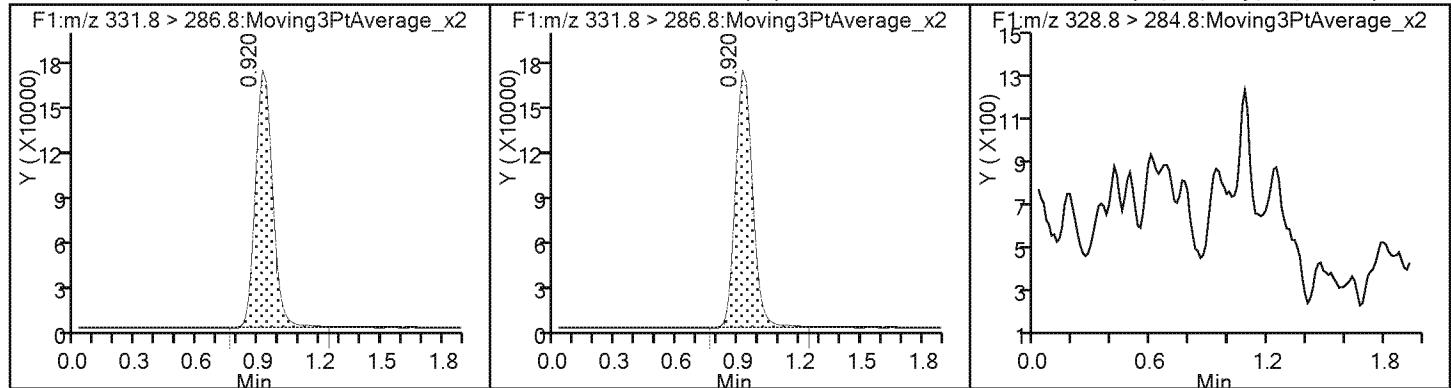

TestAmerica Denver
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30014.d
 Lims ID: 320-35355-B-6-A
 Client ID: FAY-D-701PINEA-W1-2-012418
 Sample Type: Client
 Inject. Date: 30-Jan-2018 09:48:17 ALS Bottle#: 20 Worklist Smp#: 14
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: 320-35355-B-6-A
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 30-Jan-2018 14:50:04 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK010

First Level Reviewer: meyera Date: 30-Jan-2018 14:48:31

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	12.5	125.46

FORM I
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: TestAmerica Denver Job No.: 320-35355-1

SDG No.: _____

Client Sample ID: FAY-D-663PINEA-W1-1-01241 Lab Sample ID: 320-35355-7
8

Matrix: Water Lab File ID: hfpo718A30015.d

Analysis Method: 8321A Date Collected: 01/24/2018 10:29

Extraction Method: 3535 Date Extracted: 01/29/2018 21:37

Sample wt/vol: 242.1 (mL) Date Analyzed: 01/30/2018 09:51

Con. Extract Vol.: 5 (mL) Dilution Factor: 1

Injection Volume: 20 (uL) GC Column: Synergi Hydro ID: _____

% Moisture: GPC Cleanup: (Y/N) N

Analysis Batch No.: 403246 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	0.015		0.010	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	123		50-200

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30015.d
 Lims ID: 320-35355-B-7-A
 Client ID: FAY-D-663PINEA-W1-1-012418
 Sample Type: Client
 Inject. Date: 30-Jan-2018 09:51:33 ALS Bottle#: 21 Worklist Smp#: 15
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: 320-35355-B-7-A
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 30-Jan-2018 14:50:04 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK010

First Level Reviewer: meyera Date: 30-Jan-2018 14:48:46

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 > 286.8 0.920 0.961 -0.041 1.000 902570 12.3 1448

* 2 13C3 HFPO-DA (IS)

331.8 > 286.8 0.920 0.961 -0.041 902570 10.0 1448

1 Perfluoro(2-propoxypropanoic) acid M

328.8 > 284.8 0.934 0.988 -0.054 1.000 86207 0.7278 45.6 M

QC Flag Legend

Review Flags

M - Manually Integrated

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20180130-66815.b\\hfpo718A30015.d

Injection Date: 30-Jan-2018 09:51:33

Instrument ID: LC_LCMS7

Lims ID: 320-35355-B-7-A

Lab Sample ID: 280-35355-7

Client ID: FAY-D-663PINEA-W1-1-012418

Operator ID: JBH

ALS Bottle#: 21 Worklist Smp#: 15

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A_HFPO_Du

\$ 3 13C3 HFPO-DA

* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid (M)

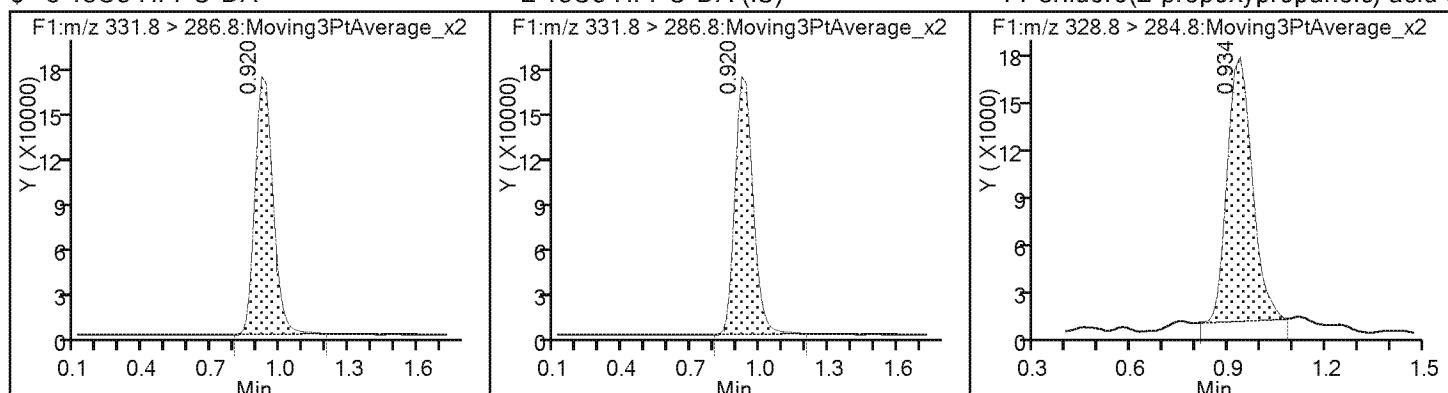

TestAmerica Denver
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30015.d
 Lims ID: 320-35355-B-7-A
 Client ID: FAY-D-663PINEA-W1-1-012418
 Sample Type: Client
 Inject. Date: 30-Jan-2018 09:51:33 ALS Bottle#: 21 Worklist Smp#: 15
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: 320-35355-B-7-A
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 30-Jan-2018 14:50:04 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM
 Process Host: XAWRK010

First Level Reviewer: meyera Date: 30-Jan-2018 14:48:46

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	12.3	123.40


TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20180130-66815.b\\hfpo718A30015.d
 Injection Date: 30-Jan-2018 09:51:33 Instrument ID: LC_LCMS7
 Lims ID: 320-35355-B-7-A Lab Sample ID: 280-35355-7
 Client ID: FAY-D-663PINEA-W1-1-012418
 Operator ID: JBH ALS Bottle#: 21 Worklist Smp#: 15
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Method: HFPO Limit Group: LC - 8321A_HFPO_Du
 Column: Detector F1:MRM

1 Perfluoro(2-propoxypropanoic) acid, CAS: 13252-13-6
 Signal: 1


RT: 0.93
 Area: 86536
 Amount: 0.731445
 Amount Units: ug/l

Processing Integration Results

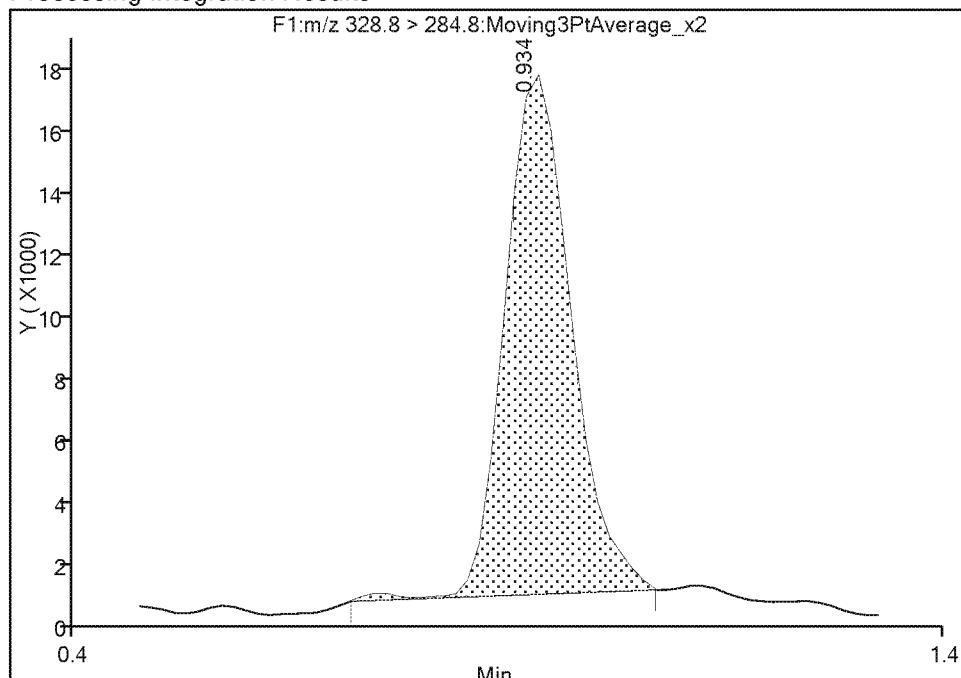

RT: 0.93
 Area: 86207
 Amount: 0.727844
 Amount Units: ug/l

Manual Integration Results

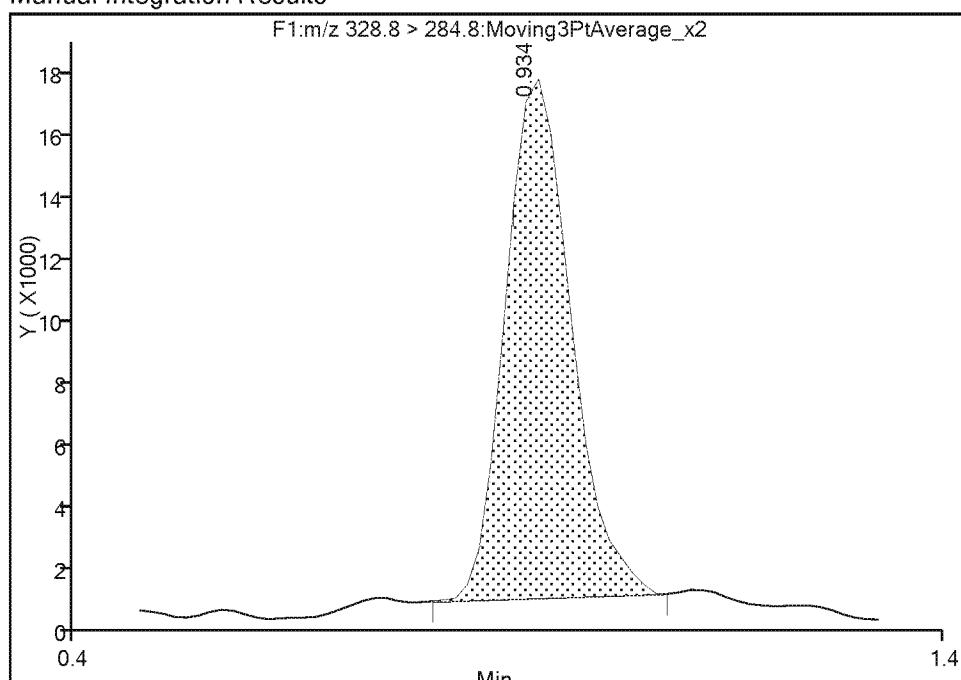

Reviewer: meyera, 30-Jan-2018 14:48:44

Audit Action: Manually Integrated

Audit Reason: Baseline

FORM I
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: TestAmerica Denver Job No.: 320-35355-1

SDG No.: _____

Client Sample ID: FAY-D-7166BUTLE-W2-1-0124 Lab Sample ID: 320-35355-8
18

Matrix: Water Lab File ID: hfpo718A30016.d

Analysis Method: 8321A Date Collected: 01/24/2018 15:43

Extraction Method: 3535 Date Extracted: 01/29/2018 21:37

Sample wt/vol: 256.1 (mL) Date Analyzed: 01/30/2018 09:54

Con. Extract Vol.: 5 (mL) Dilution Factor: 1

Injection Volume: 20 (uL) GC Column: Synergi Hydro ID: _____

% Moisture: GPC Cleanup: (Y/N) N

Analysis Batch No.: 403246 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	0.021		0.010	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	100		50-200

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30016.d
 Lims ID: 320-35355-A-8-A
 Client ID: FAY-D-7166BUTLE-W2-1-012418
 Sample Type: Client
 Inject. Date: 30-Jan-2018 09:54:48 ALS Bottle#: 22 Worklist Smp#: 16
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: 320-35355-A-8-A
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 30-Jan-2018 14:50:04 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK010


First Level Reviewer: meyera Date: 30-Jan-2018 14:48:51

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA
 331.8 > 286.8 0.920 0.961 -0.041 1.000 728977 9.97 1254
 * 2 13C3 HFPO-DA (IS)
 331.8 > 286.8 0.920 0.961 -0.041 728977 10.0 1254
 1 Perfluoro(2-propoxypropanoic) acid
 328.8 > 284.8 0.920 0.988 -0.068 1.000 94858 1.07 23.5

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20180130-66815.b\\hfpo718A30016.d
Injection Date: 30-Jan-2018 09:54:48 Instrument ID: LC_LCMS7
Lims ID: 320-35355-A-8-A Lab Sample ID: 280-35355-8
Client ID: FAY-D-7166BUTLE-W2-1-012418
Operator ID: JBH ALS Bottle#: 22 Worklist Smp#: 16
Injection Vol: 20.0 ul Dil. Factor: 1.0000
Method: HFPO Limit Group: LC - 8321A_HFPO_Du

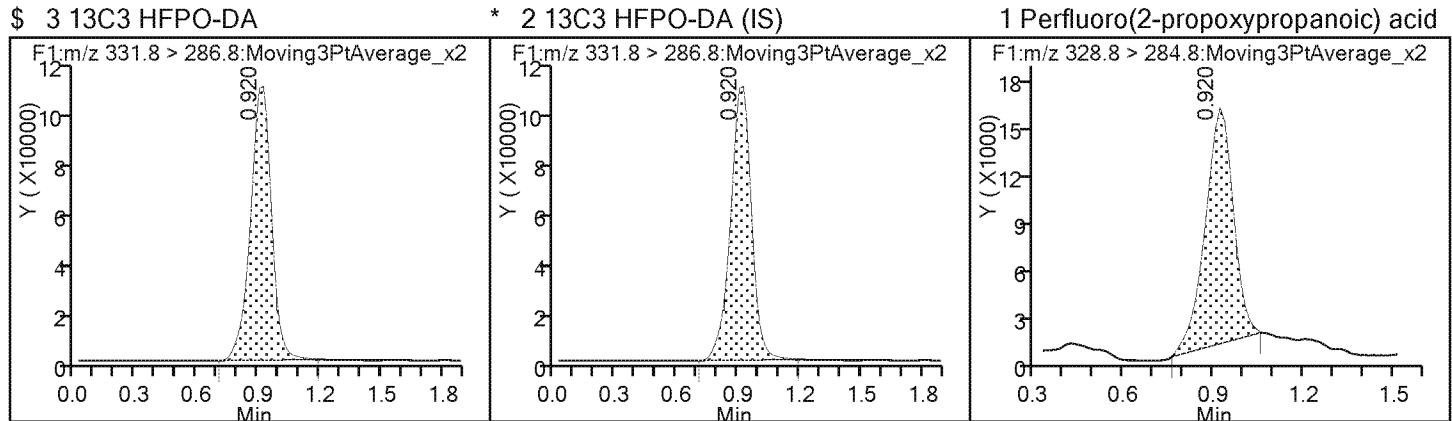

TestAmerica Denver
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30016.d
 Lims ID: 320-35355-A-8-A
 Client ID: FAY-D-7166BUTLE-W2-1-012418
 Sample Type: Client
 Inject. Date: 30-Jan-2018 09:54:48 ALS Bottle#: 22 Worklist Smp#: 16
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: 320-35355-A-8-A
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 30-Jan-2018 14:50:04 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK010

First Level Reviewer: meyera Date: 30-Jan-2018 14:48:51

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	9.97	99.66

FORM I
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: TestAmerica Denver Job No.: 320-35355-1

SDG No.: _____

Client Sample ID: FAY-D-7166BUTLE-W1-1-0124 Lab Sample ID: 320-35355-9
18

Matrix: Water Lab File ID: hfpo718A30017.d

Analysis Method: 8321A Date Collected: 01/24/2018 15:35

Extraction Method: 3535 Date Extracted: 01/29/2018 21:37

Sample wt/vol: 252.3 (mL) Date Analyzed: 01/30/2018 09:58

Con. Extract Vol.: 5 (mL) Dilution Factor: 1

Injection Volume: 20 (uL) GC Column: Synergi Hydro ID: _____

% Moisture: GPC Cleanup: (Y/N) N

Analysis Batch No.: 403246 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	0.017		0.010	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	90		50-200

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30017.d
 Lims ID: 320-35355-B-9-A
 Client ID: FAY-D-7166BUTLE-W1-1-012418
 Sample Type: Client
 Inject. Date: 30-Jan-2018 09:58:03 ALS Bottle#: 23 Worklist Smp#: 17
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: 320-35355-B-9-A
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 05-Feb-2018 13:15:57 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK025

First Level Reviewer: meyera Date: 05-Feb-2018 13:15:44

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 > 286.8 0.907 0.961 -0.054 1.000 660644 9.03 1100

* 2 13C3 HFPO-DA (IS)

331.8 > 286.8 0.907 0.961 -0.054 660644 10.0 1100

1 Perfluoro(2-propoxypropanoic) acid M

328.8 > 284.8 0.920 0.988 -0.068 1.000 72766 0.8724 11.0 M

QC Flag Legend

Review Flags

M - Manually Integrated

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20180130-66815.b\\hfpo718A30017.d

Injection Date: 30-Jan-2018 09:58:03

Instrument ID: LC_LCMS7

Lims ID: 320-35355-B-9-A

Lab Sample ID: 280-35355-9

Client ID: FAY-D-7166BUTLE-W1-1-012418

Operator ID: JBH

ALS Bottle#: 23 Worklist Smp#: 17

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A_HFPO_Du

\$ 3 13C3 HFPO-DA

* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid (M)

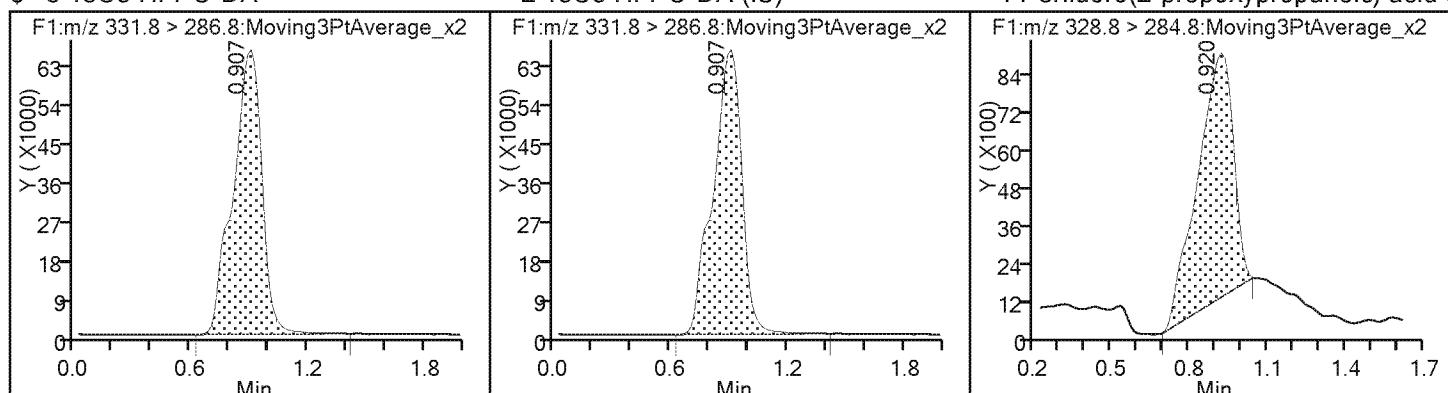

TestAmerica Denver
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30017.d
 Lims ID: 320-35355-B-9-A
 Client ID: FAY-D-7166BUTLE-W1-1-012418
 Sample Type: Client
 Inject. Date: 30-Jan-2018 09:58:03 ALS Bottle#: 23 Worklist Smp#: 17
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: 320-35355-B-9-A
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 05-Feb-2018 13:15:57 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM
 Process Host: XAWRK025

First Level Reviewer: meyera Date: 05-Feb-2018 13:15:44

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	9.03	90.32

TestAmerica Denver


Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20180130-66815.b\\hfpo718A30017.d
 Injection Date: 30-Jan-2018 09:58:03 Instrument ID: LC_LCMS7
 Lims ID: 320-35355-B-9-A Lab Sample ID: 280-35355-9
 Client ID: FAY-D-7166BUTLE-W1-1-012418
 Operator ID: JBH ALS Bottle#: 23 Worklist Smp#: 17
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Method: HFPO Limit Group: LC - 8321A_HFPO_Du
 Column: Detector F1:MRM

1 Perfluoro(2-propoxypropanoic) acid, CAS: 13252-13-6

Signal: 1


Processing Integration Results

RT: 0.92
 Area: 55586
 Amount: 0.615466
 Amount Units: ug/l

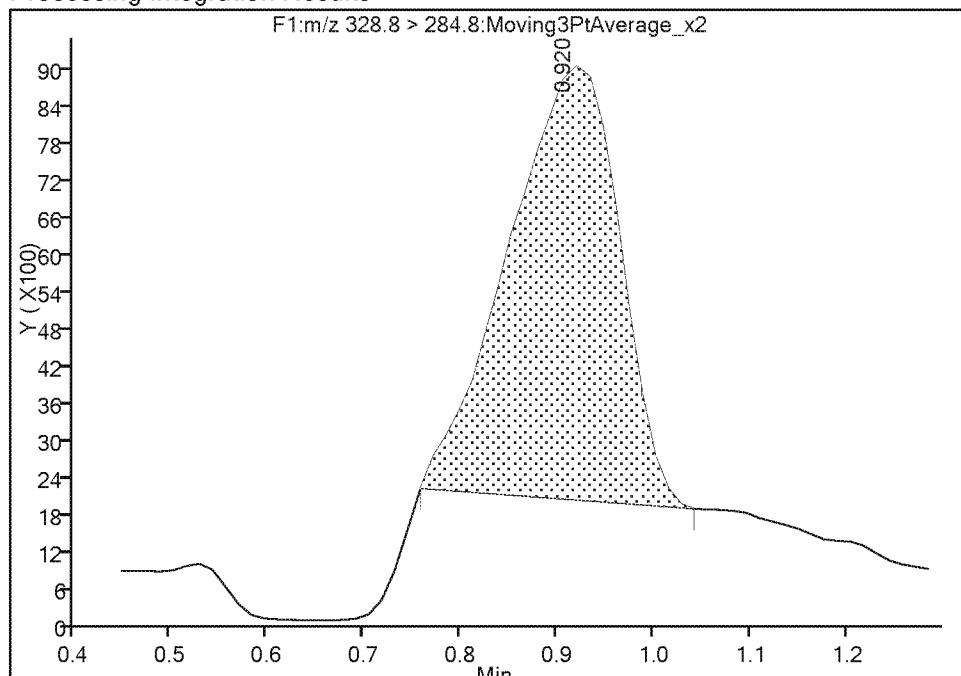

Manual Integration Results

RT: 0.92
 Area: 72766
 Amount: 0.872407
 Amount Units: ug/l

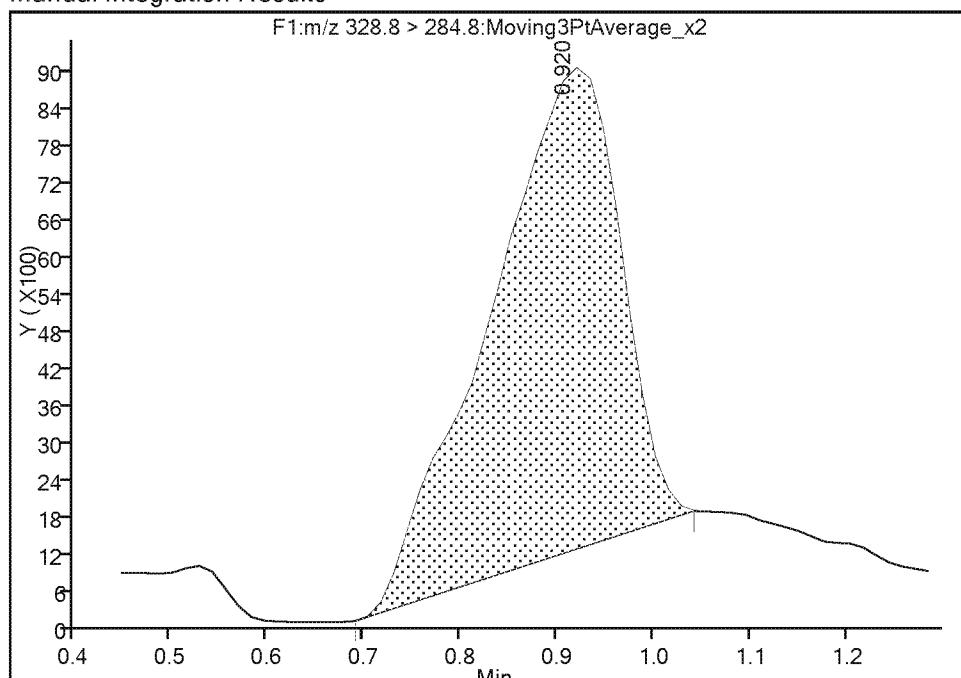

Reviewer: meyera, 05-Feb-2018 13:15:42

Audit Action: Manually Integrated

Audit Reason: Baseline

FORM I
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: TestAmerica Denver Job No.: 320-35355-1

SDG No.: _____

Client Sample ID: FAY-D-7166BUTLE-W1-1-0124 Lab Sample ID: 320-35355-10
18D _____

Matrix: Water Lab File ID: hfpo718A30020.d

Analysis Method: 8321A Date Collected: 01/24/2018 15:35

Extraction Method: 3535 Date Extracted: 01/29/2018 21:37

Sample wt/vol: 256.3 (mL) Date Analyzed: 01/30/2018 10:07

Con. Extract Vol.: 5 (mL) Dilution Factor: 1

Injection Volume: 20 (uL) GC Column: Synergi Hydro ID: _____

% Moisture: GPC Cleanup: (Y/N) N

Analysis Batch No.: 403246 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	0.021		0.010	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	103		50-200

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30020.d
 Lims ID: 320-35355-D-10-A
 Client ID: FAY-D-7166BUTLE-W1-1-012418D
 Sample Type: Client
 Inject. Date: 30-Jan-2018 10:07:48 ALS Bottle#: 26 Worklist Smp#: 20
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: 320-35355-D-10-A
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 30-Jan-2018 14:50:04 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK010

First Level Reviewer: meyera Date: 30-Jan-2018 14:49:10

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA
 331.8 > 286.8 0.907 0.961 -0.054 1.000 756680 10.3 966
 * 2 13C3 HFPO-DA (IS)
 331.8 > 286.8 0.907 0.961 -0.054 1.000 756680 10.0 966
 1 Perfluoro(2-propoxypropanoic) acid
 328.8 > 284.8 0.907 0.988 -0.081 1.000 98570 1.07 18.4


TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20180130-66815.b\\hfpo718A30020.d
Injection Date: 30-Jan-2018 10:07:48 Instrument ID: LC_LCMS7
Lims ID: 320-35355-D-10-A Lab Sample ID: 280-35355-10
Client ID: FAY-D-7166BUTLE-W1-1-012418D
Operator ID: JBH ALS Bottle#: 26 Worklist Smp#: 20
Injection Vol: 20.0 ul Dil. Factor: 1.0000
Method: HFPO Limit Group: LC - 8321A_HFPO_Du

\$ 3 13C3 HFPO-DA

* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

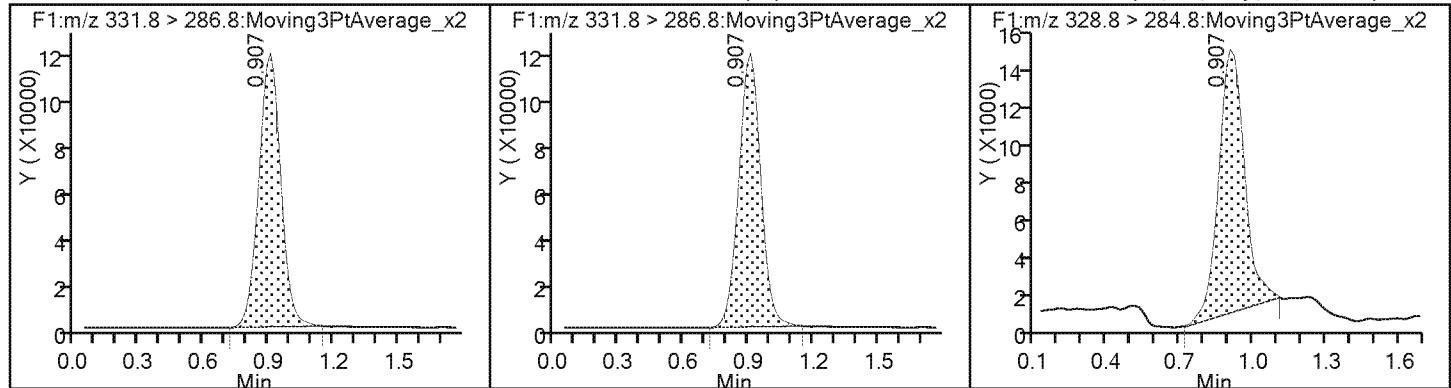

TestAmerica Denver
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30020.d
 Lims ID: 320-35355-D-10-A
 Client ID: FAY-D-7166BUTLE-W1-1-012418D
 Sample Type: Client
 Inject. Date: 30-Jan-2018 10:07:48 ALS Bottle#: 26 Worklist Smp#: 20
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: 320-35355-D-10-A
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 30-Jan-2018 14:50:04 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK010

First Level Reviewer: meyera Date: 30-Jan-2018 14:49:10

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	10.3	103.45

FORM VI
LCMS BY EXTERNAL STANDARD - INITIAL CALIBRATION DATA
RETENTION TIME SUMMARY

Lab Name: TestAmerica Denver Job No.: 320-35355-1 Analy Batch No.: 387775

SDG No.: _____

Instrument ID: LC_LCMS7 GC Column: Synergi Hyd ID: _____ Heated Purge: (Y/N) N

Calibration Start Date: 09/14/2017 14:40 Calibration End Date: 09/14/2017 15:01 Calibration ID: 30321

Calibration Files:

LEVEL:	LAB SAMPLE ID:	LAB FILE ID:
Level 1	STD001 280-387775/3	hfpo717I14052.d
Level 2	STD002 280-387775/4	hfpo717I14053.d
Level 3	STD003 280-387775/5	hfpo717I14054.d
Level 4	STD004 280-387775/6	hfpo717I14055.d
Level 5	STD005 280-387775/7	hfpo717I14056.d
Level 6	STD006 280-387775/8	hfpo717I14057.d
Level 7	STD007 280-387775/9	hfpo717I14058.d
Level 8	STD008 280-387775/10	hfpo717I14059.d

ANALYTE	LVL 1	LVL 2	LVL 3	LVL 4	LVL 5	LVL 6	LVL 7	LVL 8		RT WINDOW	AVG RT
Perfluoro(2-propoxypropanoic) acid	1.002	0.988	0.988	0.988	0.975	0.975	0.988	0.988		0.486 - 1.486	0.987
13C3 HFPO-DA	0.988	0.975	0.975	0.988	0.975	0.975	0.988	0.988		0.481 - 1.481	0.982

FORM VI
LCMS BY EXTERNAL STANDARD - INITIAL CALIBRATION DATA
CURVE EVALUATION

Lab Name: TestAmerica Denver Job No.: 320-35355-1 Analy Batch No.: 387775

SDG No.: _____

Instrument ID: LC_LCMS7 GC Column: Synergi Hyd ID: _____ Heated Purge: (Y/N) N

Calibration Start Date: 09/14/2017 14:40 Calibration End Date: 09/14/2017 15:01 Calibration ID: 30321

Calibration Files:

LEVEL:	LAB SAMPLE ID:	LAB FILE ID:
Level 1	STD001 280-387775/3	hfpo717I14052.d
Level 2	STD002 280-387775/4	hfpo717I14053.d
Level 3	STD003 280-387775/5	hfpo717I14054.d
Level 4	STD004 280-387775/6	hfpo717I14055.d
Level 5	STD005 280-387775/7	hfpo717I14056.d
Level 6	STD006 280-387775/8	hfpo717I14057.d
Level 7	STD007 280-387775/9	hfpo717I14058.d
Level 8	STD008 280-387775/10	hfpo717I14059.d

ANALYTE	CF				CURVE TYPE	COEFFICIENT			#	MIN CF	%RSD	#	MAX %RSD	R^2 OR COD	#	MIN R^2 OR COD
	LVL 1	LVL 2	LVL 3	LVL 4		B	M1	M2								
13C3 HFPO-DA	206978	200375	208177	195084	Ave		192739.525				6.4		30.0			

Note: The M1 coefficient is the same as Ave CF for an Ave curve type.

FORM VI
LCMS BY INTERNAL STANDARD - INITIAL CALIBRATION DATA
CURVE EVALUATION

Lab Name: TestAmerica Denver Job No.: 320-35355-1 Analy Batch No.: 387775

SDG No.: _____

Instrument ID: LC_LCMS7 GC Column: Synergi Hyd ID: _____ Heated Purge: (Y/N) N

Calibration Start Date: 09/14/2017 14:40 Calibration End Date: 09/14/2017 15:01 Calibration ID: 30321

ANALYTE	RRF					CURVE TYPE	COEFFICIENT			#	MIN RRF	%RSD	#	MAX %RSD	R^2 OR COD	#	MIN R^2 OR COD
	LVL 1 LVL 6	LVL 2 LVL 7	LVL 3 LVL 8	LVL 4	LVL 5		B	M1	M2								
Perfluoro(2-propoxypropanoic) acid	1.6312 0.9640	1.1780 0.9353	0.9745 0.8831	0.9868	1.0688	Lin1	0.1732	0.9076							0.9980		0.9900

Note: The M1 coefficient is the same as Ave RRF for an Ave curve type.

FORM VI
LCMS BY EXTERNAL STANDARD - INITIAL CALIBRATION DATA
RESPONSE AND CONCENTRATION

Lab Name: TestAmerica Denver Job No.: 320-35355-1 Analy Batch No.: 387775

SDG No.: _____

Instrument ID: LC_LCMS7 GC Column: Synergi Hyd ID: _____ Heated Purge: (Y/N) N

Calibration Start Date: 09/14/2017 14:40 Calibration End Date: 09/14/2017 15:01 Calibration ID: 30321

Calibration Files:

LEVEL:	LAB SAMPLE ID:	LAB FILE ID:
Level 1	STD001 280-387775/3	hfpo717I14052.d
Level 2	STD002 280-387775/4	hfpo717I14053.d
Level 3	STD003 280-387775/5	hfpo717I14054.d
Level 4	STD004 280-387775/6	hfpo717I14055.d
Level 5	STD005 280-387775/7	hfpo717I14056.d
Level 6	STD006 280-387775/8	hfpo717I14057.d
Level 7	STD007 280-387775/9	hfpo717I14058.d
Level 8	STD008 280-387775/10	hfpo717I14059.d

ANALYTE	CURVE TYPE	RESPONSE					CONCENTRATION (UG/L)				
		LVL 1 LVL 6	LVL 2 LVL 7	LVL 3 LVL 8	LVL 4	LVL 5	LVL 1 LVL 6	LVL 2 LVL 7	LVL 3 LVL 8	LVL 4	LVL 5
13C3 HFPO-DA	Ave	2069777 1724989	2003748 1884947	2081766 1878107	1950837	1824991	10.0 10.0	10.0 10.0	10.0 10.0	10.0 10.0	10.0 10.0

Curve Type Legend:

Ave = Average

FORM VI
LCMS BY INTERNAL STANDARD - INITIAL CALIBRATION DATA
RESPONSE AND CONCENTRATION

Lab Name: TestAmerica Denver Job No.: 320-35355-1 Analy Batch No.: 387775

SDG No.: _____

Instrument ID: LC_LCMS7 GC Column: Synergi Hyd ID: _____ Heated Purge: (Y/N) N

Calibration Start Date: 09/14/2017 14:40 Calibration End Date: 09/14/2017 15:01 Calibration ID: 30321

Calibration Files:

LEVEL:	LAB SAMPLE ID:	LAB FILE ID:
Level 1	STD001 280-387775/3	hfpo717I14052.d
Level 2	STD002 280-387775/4	hfpo717I14053.d
Level 3	STD003 280-387775/5	hfpo717I14054.d
Level 4	STD004 280-387775/6	hfpo717I14055.d
Level 5	STD005 280-387775/7	hfpo717I14056.d
Level 6	STD006 280-387775/8	hfpo717I14057.d
Level 7	STD007 280-387775/9	hfpo717I14058.d
Level 8	STD008 280-387775/10	hfpo717I14059.d

ANALYTE	IS REF	CURVE TYPE	RESPONSE					CONCENTRATION (UG/L)				
			LVL 1 LVL 6	LVL 2 LVL 7	LVL 3 LVL 8	LVL 4	LVL 5	LVL 1 LVL 6	LVL 2 LVL 7	LVL 3 LVL 8	LVL 4	LVL 5
Perfluoro(2-propoxypropanoic) acid	13CP ODA	Lin1	84406 1662919	118017 4407541	202876 8293101	385009	975278	0.250 10.0	0.500 25.0	1.00 50.0	2.00	5.00

Curve Type Legend:

Lin1 = Linear 1/conc ISTD

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20170915-62647.b\hfpo717I14052.d
 Lims ID: std001
 Client ID:
 Sample Type: IC Calib Level: 1
 Inject. Date: 14-Sep-2017 14:40:03 ALS Bottle#: 2 Worklist Smp#: 3
 Injection Vol: 10.0 ul Dil. Factor: 1.0000
 Sample Info: L1
 Misc. Info.: HFPO17I14
 Operator ID: JBH Instrument ID: LC_LCMS7
 Sublist: chrom-HFPO*sub1

Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20170915-62647.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 15-Sep-2017 07:29:39 Calib Date: 14-Sep-2017 15:01:22
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20170915-62647.b\hfpo717I14059.d

Column 1 : Det: F1:MRM
 Process Host: XAWRK034

First Level Reviewer: meyera Date: 15-Sep-2017 07:28:15

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA
 331.8 > 286.8 0.988 0.981 0.007 1.000 2069777 10.7 429
 * 2 13C3 HFPO-DA (IS)
 331.8 > 286.8 0.988 0.981 0.007 2069777 10.0 429
 1 Perfluoro(2-propoxypropanoic) acid
 328.8 > 284.8 1.002 0.986 0.016 1.000 84406 0.2585 49.7

Reagents:

HFPO_CAL-1_00030 Amount Added: 1.00 Units: mL

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20170915-62647.b\\hfp0717\\14052.d

Injection Date: 14-Sep-2017 14:40:03 Instrument ID: LC_LCMS7

Lims ID: std001

Client ID:

Operator ID: JBH ALS Bottle#: 2 Worklist Smp#: 3


Injection Vol: 10.0 ul Dil. Factor: 1.0000

Method: HFPO Limit Group: LC - 8321A_HFPO_Du

\$ 3 13C3 HFPO-DA

* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

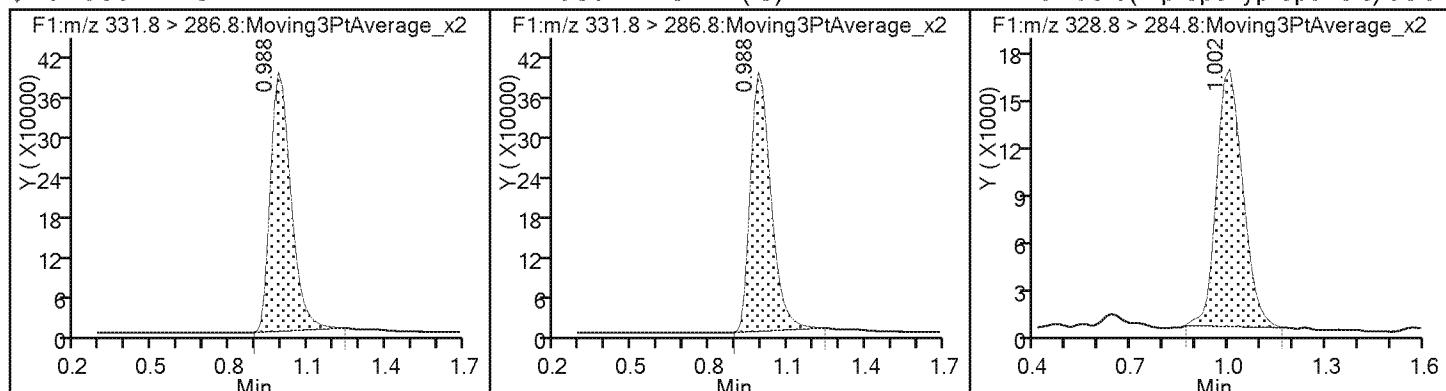

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20170915-62647.b\hfpo717I14053.d
 Lims ID: std002
 Client ID:
 Sample Type: IC Calib Level: 2
 Inject. Date: 14-Sep-2017 14:43:06 ALS Bottle#: 3 Worklist Smp#: 4
 Injection Vol: 10.0 ul Dil. Factor: 1.0000
 Sample Info: L2
 Misc. Info.: HFPO17I14
 Operator ID: JBH Instrument ID: LC_LCMS7
 Sublist: chrom-HFPO*sub1

Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20170915-62647.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 15-Sep-2017 07:29:39 Calib Date: 14-Sep-2017 15:01:22
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20170915-62647.b\hfpo717I14059.d

Column 1 : Det: F1:MRM

Process Host: XAWRK034

First Level Reviewer: meyera Date: 15-Sep-2017 07:28:18

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

* 2 13C3 HFPO-DA (IS)

331.8 > 286.8 0.975 0.981 -0.006 2003748 10.0 386

\$ 3 13C3 HFPO-DA

331.8 > 286.8 0.975 0.981 -0.006 1.000 2003748 10.4 386

1 Perfluoro(2-propoxypropanoic) acid

328.8 > 284.8 0.988 0.986 0.002 1.000 118017 0.4581 56.6

Reagents:

HFPO_CAL-2_00031 Amount Added: 1.00 Units: mL

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20170915-62647.b\\hfp0717\\14053.d

Injection Date: 14-Sep-2017 14:43:06

Instrument ID: LC_LCMS7

Lims ID: std002

Client ID:

Operator ID: JBH

ALS Bottle#: 3 Worklist Smp#: 4

Injection Vol: 10.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A_HFPO_Du

* 2 13C3 HFPO-DA (IS)

\$ 3 13C3 HFPO-DA

1 Perfluoro(2-propoxypropanoic) acid

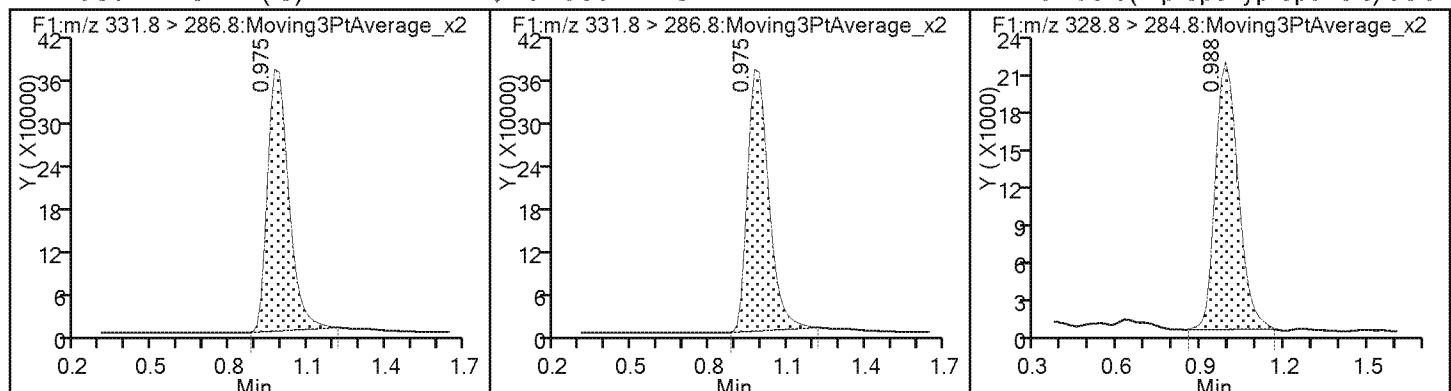

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20170915-62647.b\hfpo717I14054.d
 Lims ID: std003
 Client ID:
 Sample Type: IC Calib Level: 3
 Inject. Date: 14-Sep-2017 14:46:08 ALS Bottle#: 4 Worklist Smp#: 5
 Injection Vol: 10.0 ul Dil. Factor: 1.0000
 Sample Info: L3
 Misc. Info.: HFPO17I14
 Operator ID: JBH Instrument ID: LC_LCMS7
 Sublist: chrom-HFPO*sub1

Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20170915-62647.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 15-Sep-2017 07:29:40 Calib Date: 14-Sep-2017 15:01:22
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20170915-62647.b\hfpo717I14059.d

Column 1 : Det: F1:MRM

Process Host: XAWRK034

First Level Reviewer: meyera Date: 15-Sep-2017 07:28:20

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA
331.8 > 286.8 0.975 0.981 -0.006 1.000 2081766 10.8 403

* 2 13C3 HFPO-DA (IS)
331.8 > 286.8 0.975 0.981 -0.006 2081766 10.0 403

1 Perfluoro(2-propoxypropanoic) acid
328.8 > 284.8 0.988 0.986 0.002 1.000 202876 0.8830 108

Reagents:

HFPO_CAL-3_00030 Amount Added: 1.00 Units: mL

Report Date: 15-Sep-2017 07:29:40

Chrom Revision: 2.2 16-Aug-2017 16:24:46

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20170915-62647.b\\hfp0717\\14054.d

Injection Date: 14-Sep-2017 14:46:08

Instrument ID: LC_LCMS7

Lims ID: std003

Client ID:

Operator ID: JBH

ALS Bottle#: 4 Worklist Smp#: 5

Injection Vol: 10.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A_HFPO_Du

\$ 3 13C3 HFPO-DA

* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

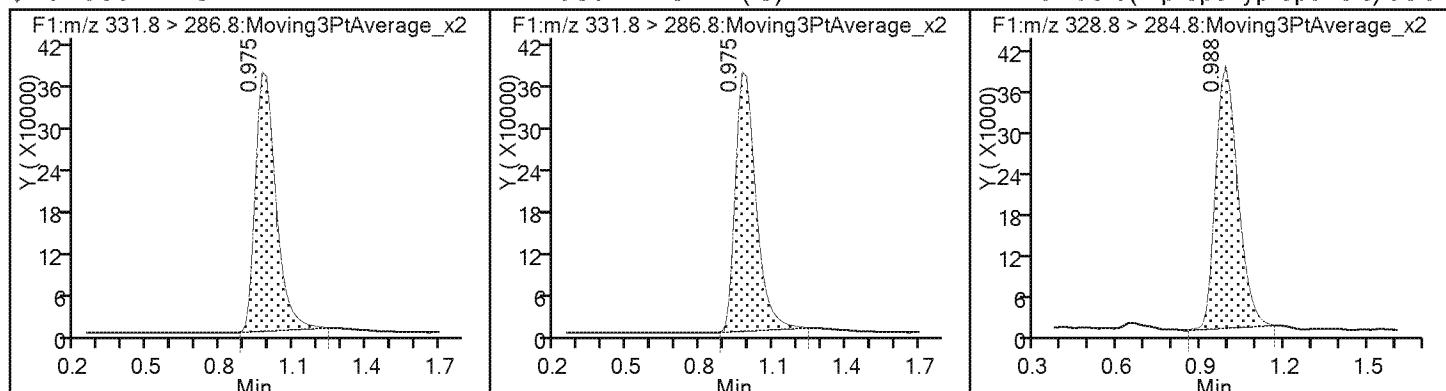

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20170915-62647.b\hfpo717I14055.d
 Lims ID: std004
 Client ID:
 Sample Type: IC Calib Level: 4
 Inject. Date: 14-Sep-2017 14:49:11 ALS Bottle#: 5 Worklist Smp#: 6
 Injection Vol: 10.0 ul Dil. Factor: 1.0000
 Sample Info: L4
 Misc. Info.: HFPO17I14
 Operator ID: JBH Instrument ID: LC_LCMS7
 Sublist: chrom-HFPO*sub1

Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20170915-62647.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 15-Sep-2017 07:29:40 Calib Date: 14-Sep-2017 15:01:22
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20170915-62647.b\hfpo717I14059.d

Column 1 : Det: F1:MRM
 Process Host: XAWRK034

First Level Reviewer: meyera Date: 15-Sep-2017 07:28:22

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

* 2 13C3 HFPO-DA (IS)
 331.8 > 286.8 0.988 0.981 0.007 1950837 10.0 384
 \$ 3 13C3 HFPO-DA
 331.8 > 286.8 0.988 0.981 0.007 1.000 1950837 10.1 384
 1 Perfluoro(2-propoxypropanoic) acid
 328.8 > 284.8 0.988 0.986 0.002 1.000 385009 1.98 162

Reagents:

HFPO_CAL-4_00030 Amount Added: 1.00 Units: mL

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20170915-62647.b\\hfp0717\\14055.d

Injection Date: 14-Sep-2017 14:49:11 Instrument ID: LC_LCMS7

Lims ID: std004

Client ID:

Operator ID: JBH ALS Bottle#: 5 Worklist Smp#: 6


Injection Vol: 10.0 ul Dil. Factor: 1.0000

Method: HFPO Limit Group: LC - 8321A_HFPO_Du

* 2 13C3 HFPO-DA (IS)

\$ 3 13C3 HFPO-DA

1 Perfluoro(2-propoxypropanoic) acid

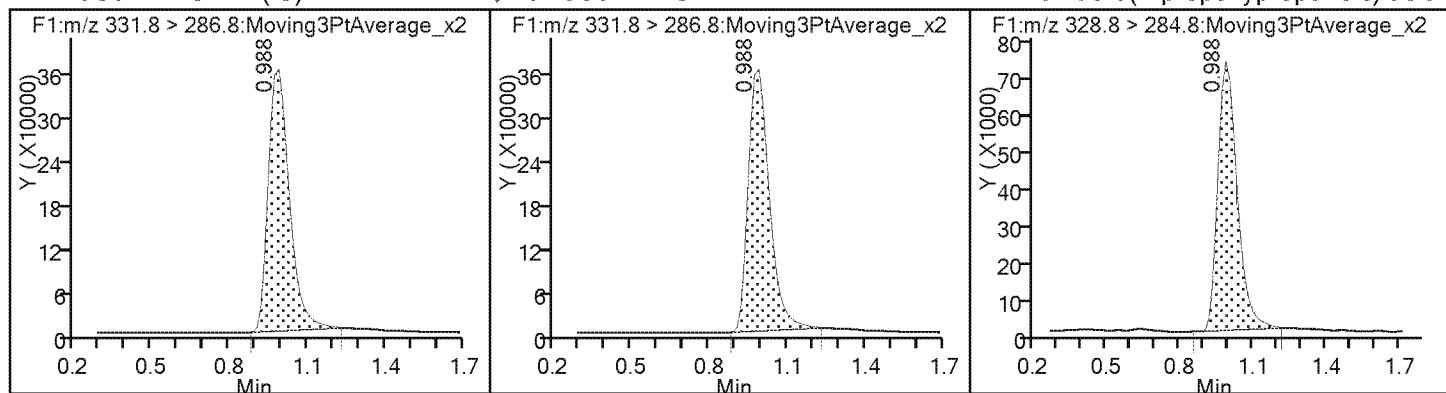

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20170915-62647.b\hfpo717I14056.d
 Lims ID: std005
 Client ID:
 Sample Type: IC Calib Level: 5
 Inject. Date: 14-Sep-2017 14:52:13 ALS Bottle#: 6 Worklist Smp#: 7
 Injection Vol: 10.0 ul Dil. Factor: 1.0000
 Sample Info: L5
 Misc. Info.: HFPO17I14
 Operator ID: JBH Instrument ID: LC_LCMS7
 Sublist: chrom-HFPO*sub1

Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20170915-62647.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 15-Sep-2017 07:29:41 Calib Date: 14-Sep-2017 15:01:22
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20170915-62647.b\hfpo717I14059.d

Column 1 : Det: F1:MRM
 Process Host: XAWRK034

First Level Reviewer: meyera Date: 15-Sep-2017 07:28:25

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA
 331.8 > 286.8 0.975 0.981 -0.006 1.000 1824991 9.47 371
 * 2 13C3 HFPO-DA (IS)
 331.8 > 286.8 0.975 0.981 -0.006 1824991 10.0 371
 1 Perfluoro(2-propoxypropanoic) acid
 328.8 > 284.8 0.975 0.986 -0.011 1.000 975278 5.70 268

Reagents:

HFPO_CAL-5_00067 Amount Added: 1.00 Units: mL

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20170915-62647.b\\hfp0717\\14056.d

Injection Date: 14-Sep-2017 14:52:13 Instrument ID: LC_LCMS7

Lims ID: std005

Client ID:

Operator ID: JBH ALS Bottle#: 6 Worklist Smp#: 7


Injection Vol: 10.0 ul Dil. Factor: 1.0000

Method: HFPO Limit Group: LC - 8321A_HFPO_Du

\$ 3 13C3 HFPO-DA

* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

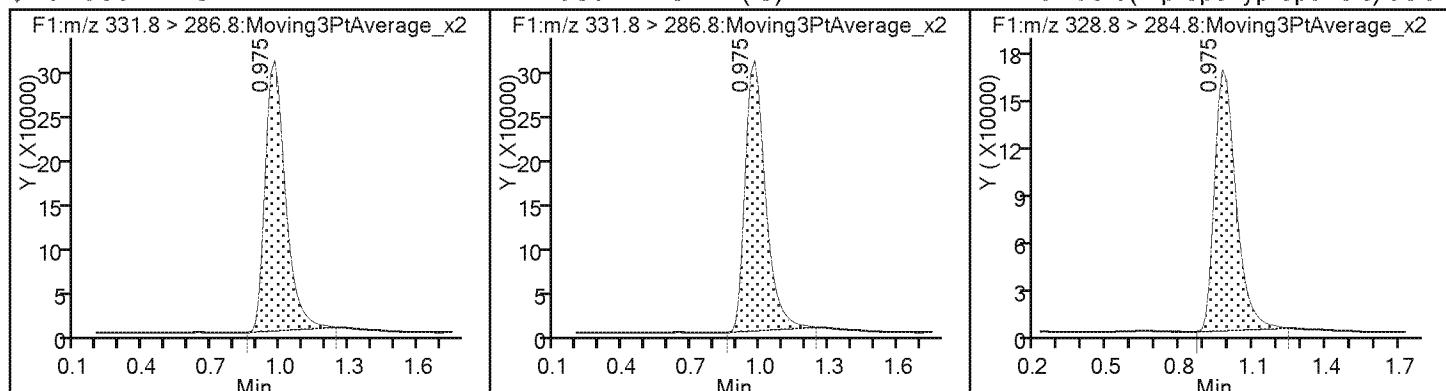

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20170915-62647.b\hfpo717I14057.d
 Lims ID: std006
 Client ID:
 Sample Type: IC Calib Level: 6
 Inject. Date: 14-Sep-2017 14:55:16 ALS Bottle#: 7 Worklist Smp#: 8
 Injection Vol: 10.0 ul Dil. Factor: 1.0000
 Sample Info: L6
 Misc. Info.: HFPO17I14
 Operator ID: JBH Instrument ID: LC_LCMS7
 Sublist: chrom-HFPO*sub1

Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20170915-62647.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 15-Sep-2017 07:29:41 Calib Date: 14-Sep-2017 15:01:22
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20170915-62647.b\hfpo717I14059.d

Column 1 : Det: F1:MRM

Process Host: XAWRK034

First Level Reviewer: meyera Date: 15-Sep-2017 07:28:27

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

* 2 13C3 HFPO-DA (IS)

331.8 > 286.8 0.975 0.981 -0.006 1724989 10.0 287

\$ 3 13C3 HFPO-DA

331.8 > 286.8 0.975 0.981 -0.006 1.000 1724989 8.95 287

1 Perfluoro(2-propoxypropanoic) acid

328.8 > 284.8 0.975 0.986 -0.011 1.000 1662919 10.4 248

Reagents:

HFPO_CAL-6_00067 Amount Added: 1.00 Units: mL

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20170915-62647.b\\hfp0717\\14057.d

Injection Date: 14-Sep-2017 14:55:16 Instrument ID: LC_LCMS7

Lims ID: std006

Client ID:

Operator ID: JBH ALS Bottle#: 7 Worklist Smp#: 8


Injection Vol: 10.0 ul Dil. Factor: 1.0000

Method: HFPO Limit Group: LC - 8321A_HFPO_Du

* 2 13C3 HFPO-DA (IS)

\$ 3 13C3 HFPO-DA

1 Perfluoro(2-propoxypropanoic) acid

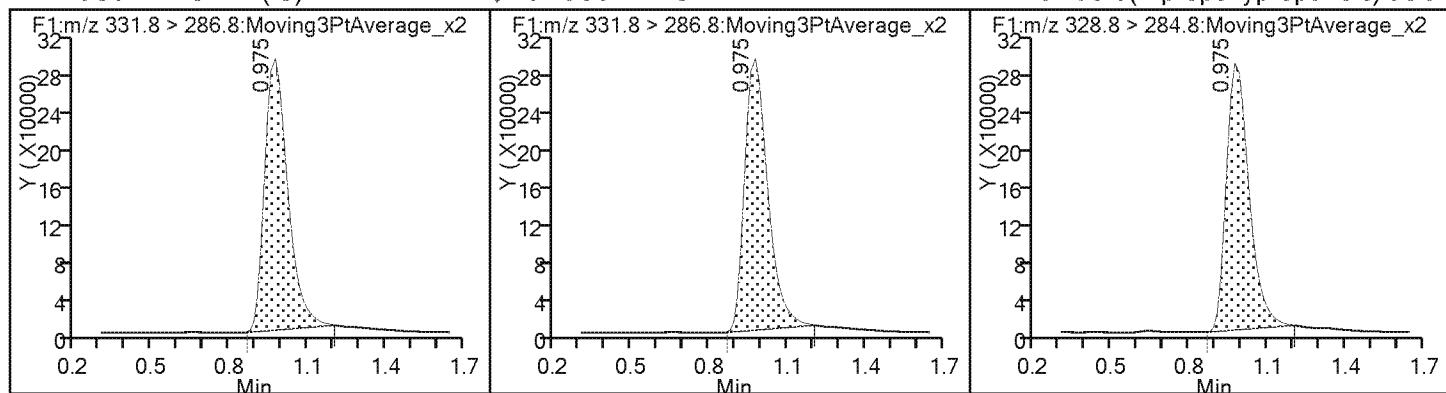

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20170915-62647.b\hfpo717I14058.d
 Lims ID: std007
 Client ID:
 Sample Type: IC Calib Level: 7
 Inject. Date: 14-Sep-2017 14:58:19 ALS Bottle#: 8 Worklist Smp#: 9
 Injection Vol: 10.0 ul Dil. Factor: 1.0000
 Sample Info: L7
 Misc. Info.: HFPO17I14
 Operator ID: JBH Instrument ID: LC_LCMS7
 Sublist: chrom-HFPO*sub1

Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20170915-62647.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 15-Sep-2017 07:29:42 Calib Date: 14-Sep-2017 15:01:22
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20170915-62647.b\hfpo717I14059.d

Column 1 : Det: F1:MRM
 Process Host: XAWRK034

First Level Reviewer: meyera Date: 15-Sep-2017 07:28:29

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA
 331.8 > 286.8 0.988 0.981 0.007 1.000 1884947 9.78 361
 * 2 13C3 HFPO-DA (IS)
 331.8 > 286.8 0.988 0.981 0.007 1884947 10.0 361
 1 Perfluoro(2-propoxypropanoic) acid
 328.8 > 284.8 0.988 0.986 0.002 1.000 4407541 25.6 379

Reagents:

HFPO_CAL-7_00030 Amount Added: 1.00 Units: mL

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20170915-62647.b\\hfpo717\\14058.d

Injection Date: 14-Sep-2017 14:58:19 Instrument ID: LC_LCMS7

Lims ID: std007

Client ID:

Operator ID: JBH ALS Bottle#: 8 Worklist Smp#: 9


Injection Vol: 10.0 ul Dil. Factor: 1.0000

Method: HFPO Limit Group: LC - 8321A_HFPO_Du

\$ 3 13C3 HFPO-DA

* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

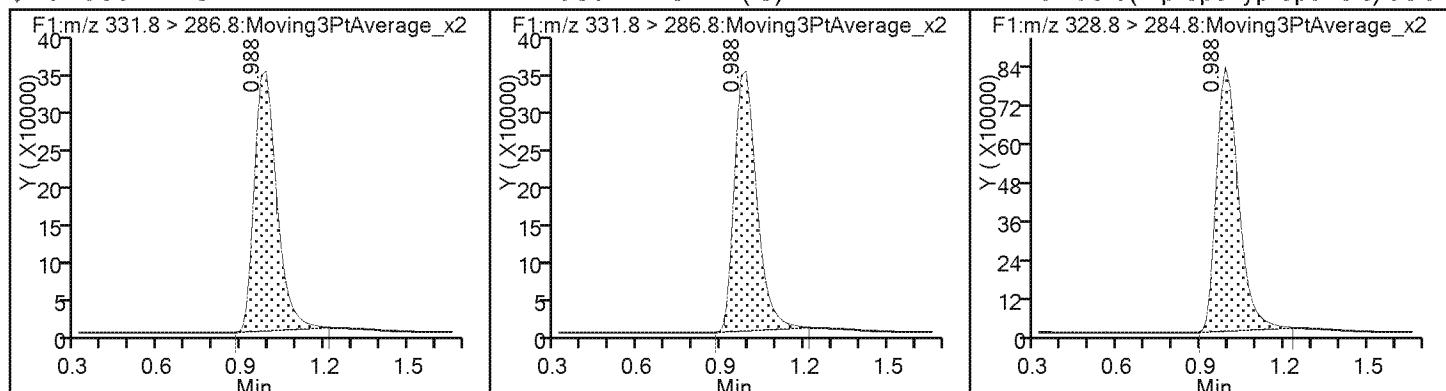

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20170915-62647.b\hfpo717I14059.d
 Lims ID: std008
 Client ID:
 Sample Type: IC Calib Level: 8
 Inject. Date: 14-Sep-2017 15:01:22 ALS Bottle#: 9 Worklist Smp#: 10
 Injection Vol: 10.0 ul Dil. Factor: 1.0000
 Sample Info: L8
 Misc. Info.: HFPO17I14
 Operator ID: JBH Instrument ID: LC_LCMS7
 Sublist: chrom-HFPO*sub1

Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20170915-62647.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 15-Sep-2017 07:29:43 Calib Date: 14-Sep-2017 15:01:22
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20170915-62647.b\hfpo717I14059.d

Column 1 : Det: F1:MRM
 Process Host: XAWRK034

First Level Reviewer: meyera Date: 15-Sep-2017 07:28:32

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

* 2 13C3 HFPO-DA (IS)
 331.8 > 286.8 0.988 0.981 0.007 1878107 10.0 379

\$ 3 13C3 HFPO-DA
 331.8 > 286.8 0.988 0.981 0.007 1.000 1878107 9.74 379

1 Perfluoro(2-propoxypropanoic) acid
 328.8 > 284.8 0.988 0.986 0.002 1.000 8293101 48.5 359

Reagents:

HFPO_CAL-8_00030 Amount Added: 1.00 Units: mL

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20170915-62647.b\\hfpo717\\14059.d

Injection Date: 14-Sep-2017 15:01:22 Instrument ID: LC_LCMS7

Lims ID: std008

Client ID:

Operator ID: JBH

ALS Bottle#: 9 Worklist Smp#: 10

Injection Vol: 10.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A_HFPO_Du

* 2 13C3 HFPO-DA (IS)

\$ 3 13C3 HFPO-DA

1 Perfluoro(2-propoxypropanoic) acid

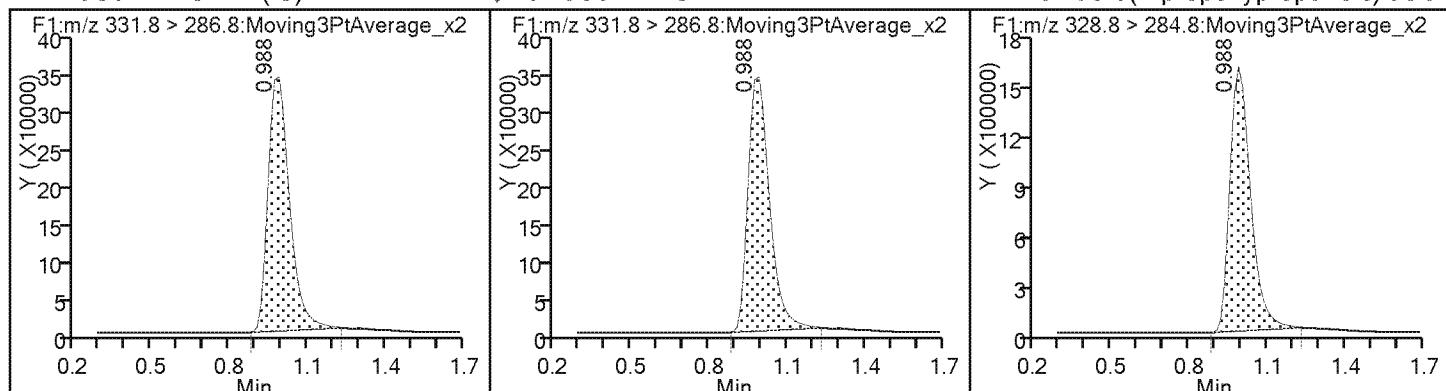

FORM VI
LCMS BY EXTERNAL STANDARD - INITIAL CALIBRATION DATA
RETENTION TIME SUMMARY

Lab Name: TestAmerica Denver Job No.: 320-35355-1 Analy Batch No.: 390728
SDG No.: _____
Instrument ID: LC_LCMS7 GC Column: Synergi Hyd ID: _____ Heated Purge: (Y/N) N
Calibration Start Date: 10/10/2017 09:35 Calibration End Date: 10/10/2017 09:58 Calibration ID: 30558

Calibration Files:

LEVEL:	LAB SAMPLE ID:	LAB FILE ID:
Level 1	STD001 280-390728/3	hfpo717J10026.d
Level 2	STD002 280-390728/4	hfpo717J10027.d
Level 3	STD003 280-390728/5	hfpo717J10028.d
Level 4	STD004 280-390728/6	hfpo717J10029.d
Level 5	STD005 280-390728/7	hfpo717J10030.d
Level 6	STD006 280-390728/8	hfpo717J10031.d
Level 7	STD007 280-390728/9	hfpo717J10032.d
Level 8	STD008 280-390728/10	hfpo717J10033.d

ANALYTE	LVL 1	LVL 2	LVL 3	LVL 4	LVL 5	LVL 6	LVL 7	LVL 8		RT WINDOW	AVG RT
HFPO-DA	0.893	0.880	0.880	0.880	0.893	0.880	0.880	0.893		0.385 - 1.385	0.885
13C3 HFPO-DA	0.880	0.880	0.880	0.880	0.880	0.880	0.880	0.880		0.380 - 1.380	0.880

FORM VI
LCMS BY EXTERNAL STANDARD - INITIAL CALIBRATION DATA
CURVE EVALUATION

Lab Name: TestAmerica Denver

Job No.: 320-35355-1

Analy Batch No.: 390728

SDG No.:

Instrument ID: LC_LCMS7 GC Column: Synergi Hyd ID: Heated Purge: (Y/N) N

Calibration Start Date: 10/10/2017 09:35 Calibration End Date: 10/10/2017 09:58 Calibration ID: 30558

Calibration Files:

LEVEL:	LAB SAMPLE ID:	LAB FILE ID:
Level 1	STD001 280-390728/3	hfpo717J10026.d
Level 2	STD002 280-390728/4	hfpo717J10027.d
Level 3	STD003 280-390728/5	hfpo717J10028.d
Level 4	STD004 280-390728/6	hfpo717J10029.d
Level 5	STD005 280-390728/7	hfpo717J10030.d
Level 6	STD006 280-390728/8	hfpo717J10031.d
Level 7	STD007 280-390728/9	hfpo717J10032.d
Level 8	STD008 280-390728/10	hfpo717J10033.d

ANALYTE	CF				CURVE TYPE	COEFFICIENT			#	MIN CF	%RSD	#	MAX %RSD	R^2 OR COD	#	MIN R^2 OR COD
	LVL 1	LVL 2	LVL 3	LVL 4		B	M1	M2								
13C3 HFPO-DA	73075 74460	74523 73194	75043 72919	71803 70142	Ave		73144.6750				2.2	30.0				

Note: The M1 coefficient is the same as Ave CF for an Ave curve type.

FORM VI
LCMS BY INTERNAL STANDARD - INITIAL CALIBRATION DATA
CURVE EVALUATION

Lab Name: TestAmerica Denver Job No.: 320-35355-1 Analy Batch No.: 390728

SDG No.: _____

Instrument ID: LC_LCMS7 GC Column: Synergi Hyd ID: _____ Heated Purge: (Y/N) N

Calibration Start Date: 10/10/2017 09:35 Calibration End Date: 10/10/2017 09:58 Calibration ID: 30558

ANALYTE	RRF					CURVE TYPE	COEFFICIENT			#	MIN RRF	%RSD	#	MAX %RSD	R^2 OR COD	#	MIN R^2 OR COD
	LVL 1 LVL 6	LVL 2 LVL 7	LVL 3 LVL 8	LVL 4	LVL 5		B	M1	M2								
HFPO-DA	1.6980 1.0102	1.7128 0.9824	1.1896 1.0419	1.1637	1.0154	Lin1	0.2185	1.0121							0.9980		0.9900

Note: The M1 coefficient is the same as Ave RRF for an Ave curve type.

FORM VI
LCMS BY EXTERNAL STANDARD - INITIAL CALIBRATION DATA
RESPONSE AND CONCENTRATION

Lab Name: TestAmerica Denver Job No.: 320-35355-1 Analy Batch No.: 390728
SDG No.: _____
Instrument ID: LC_LCMS7 GC Column: Synergi Hyd ID: _____ Heated Purge: (Y/N) N
Calibration Start Date: 10/10/2017 09:35 Calibration End Date: 10/10/2017 09:58 Calibration ID: 30558

Calibration Files:

LEVEL:	LAB SAMPLE ID:	LAB FILE ID:
Level 1	STD001 280-390728/3	hfpo717J10026.d
Level 2	STD002 280-390728/4	hfpo717J10027.d
Level 3	STD003 280-390728/5	hfpo717J10028.d
Level 4	STD004 280-390728/6	hfpo717J10029.d
Level 5	STD005 280-390728/7	hfpo717J10030.d
Level 6	STD006 280-390728/8	hfpo717J10031.d
Level 7	STD007 280-390728/9	hfpo717J10032.d
Level 8	STD008 280-390728/10	hfpo717J10033.d

ANALYTE	CURVE TYPE	RESPONSE					CONCENTRATION (UG/L)				
		LVL 1 LVL 6	LVL 2 LVL 7	LVL 3 LVL 8	LVL 4	LVL 5	LVL 1 LVL 6	LVL 2 LVL 7	LVL 3 LVL 8	LVL 4	LVL 5
13C3 HFPO-DA	Ave	730749 731935	745227 729188	750427 701420	718028	744600	10.0 10.0	10.0 10.0	10.0 10.0	10.0	10.0

Curve Type Legend:

Ave = Average

FORM VI
LCMS BY INTERNAL STANDARD - INITIAL CALIBRATION DATA
RESPONSE AND CONCENTRATION

Lab Name: TestAmerica Denver Job No.: 320-35355-1 Analy Batch No.: 390728

SDG No.: _____

Instrument ID: LC_LCMS7 GC Column: Synergi Hyd ID: _____ Heated Purge: (Y/N) N

Calibration Start Date: 10/10/2017 09:35 Calibration End Date: 10/10/2017 09:58 Calibration ID: 30558

Calibration Files:

LEVEL:	LAB SAMPLE ID:	LAB FILE ID:
Level 1	STD001 280-390728/3	hfpo717J10026.d
Level 2	STD002 280-390728/4	hfpo717J10027.d
Level 3	STD003 280-390728/5	hfpo717J10028.d
Level 4	STD004 280-390728/6	hfpo717J10029.d
Level 5	STD005 280-390728/7	hfpo717J10030.d
Level 6	STD006 280-390728/8	hfpo717J10031.d
Level 7	STD007 280-390728/9	hfpo717J10032.d
Level 8	STD008 280-390728/10	hfpo717J10033.d

ANALYTE	IS REF	CURVE TYPE	RESPONSE					CONCENTRATION (UG/L)				
			LVL 1 LVL 6	LVL 2 LVL 7	LVL 3 LVL 8	LVL 4	LVL 5	LVL 1 LVL 6	LVL 2 LVL 7	LVL 3 LVL 8	LVL 4	LVL 5
HFPO-DA	13CP ODA	Lin1	31020 739399	63823 1790812	89272 3654104	167109	378047	0.250 10.0	0.500 25.0	1.00 50.0	2.00	5.00

Curve Type Legend:

Lin1 = Linear 1/conc ISTD

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10026.d
 Lims ID: std001
 Client ID:
 Sample Type: IC Calib Level: 1
 Inject. Date: 10-Oct-2017 09:35:28 ALS Bottle#: 2 Worklist Smp#: 3
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: L1
 Misc. Info.: HFPO17J10
 Operator ID: JBH Instrument ID: LC_LCMS7
 Sublist: chrom-HFPO*sub1

Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 10-Oct-2017 12:51:45 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM
 Process Host: XAWRK005

First Level Reviewer: meyera Date: 10-Oct-2017 11:50:42

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA
 331.8 > 286.8 0.880 0.880 0.0 1.000 730749 10.0 397
 * 2 13C3 HFPO-DA (IS)
 331.8 > 286.8 0.880 0.880 0.0 1.000 730749 10.0 397
 1 Perfluoro(2-propoxypropanoic) acid M
 328.8 > 284.8 0.893 0.885 0.008 1.000 31020 0.2036 14.1 M

QC Flag Legend

Review Flags

M - Manually Integrated

Reagents:

HFPO_CAL-1_00031 Amount Added: 1.00 Units: mL

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20171010-63483.b\\hfp0717J10026.d

Injection Date: 10-Oct-2017 09:35:28

Instrument ID: LC_LCMS7

Lims ID: std001

Client ID:

Operator ID: JBH

ALS Bottle#: 2 Worklist Smp#: 3

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A_HFPO_Du

\$ 3 13C3 HFPO-DA

* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid (M)

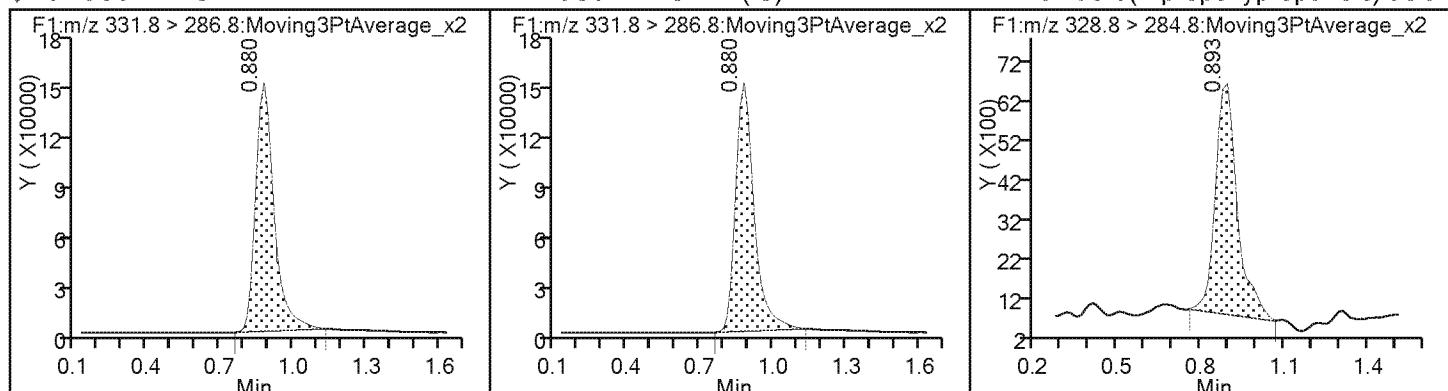

TestAmerica Denver


Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20171010-63483.b\\hfpo717J10026.d
 Injection Date: 10-Oct-2017 09:35:28 Instrument ID: LC_LCMS7
 Lims ID: std001
 Client ID:
 Operator ID: JBH ALS Bottle#: 2 Worklist Smp#: 3
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Method: HFPO Limit Group: LC - 8321A_HFPO_Du
 Column: Detector F1:MRM

1 Perfluoro(2-propoxypropanoic) acid, CAS: 13252-13-6

Signal: 1


RT: 0.89
 Area: 24407
 Amount: 0.162386
 Amount Units: ug/l

Processing Integration Results

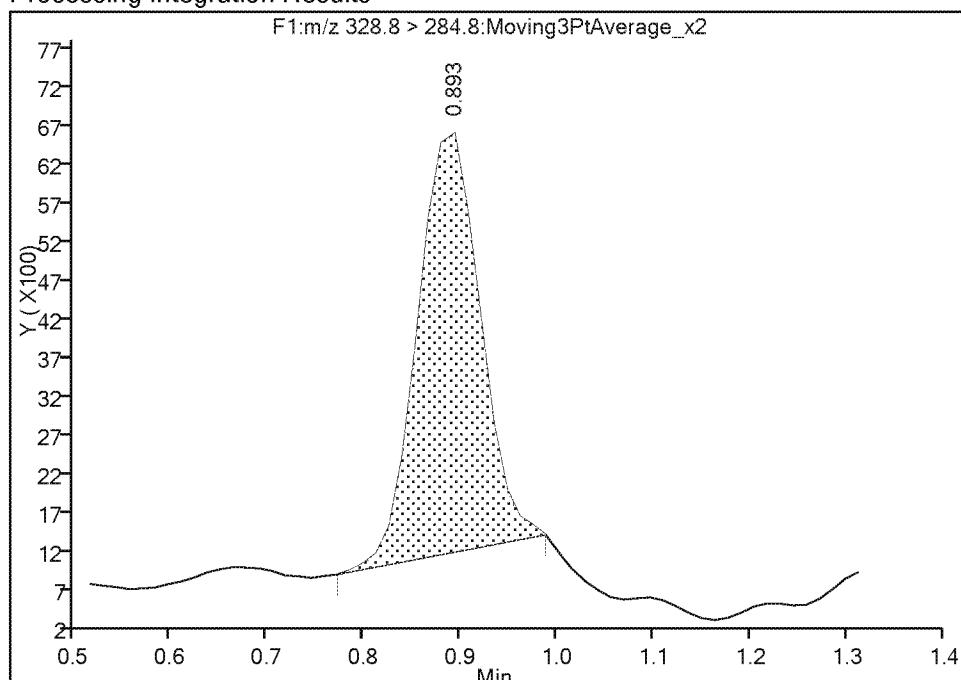

RT: 0.89
 Area: 31020
 Amount: 0.203553
 Amount Units: ug/l

Manual Integration Results

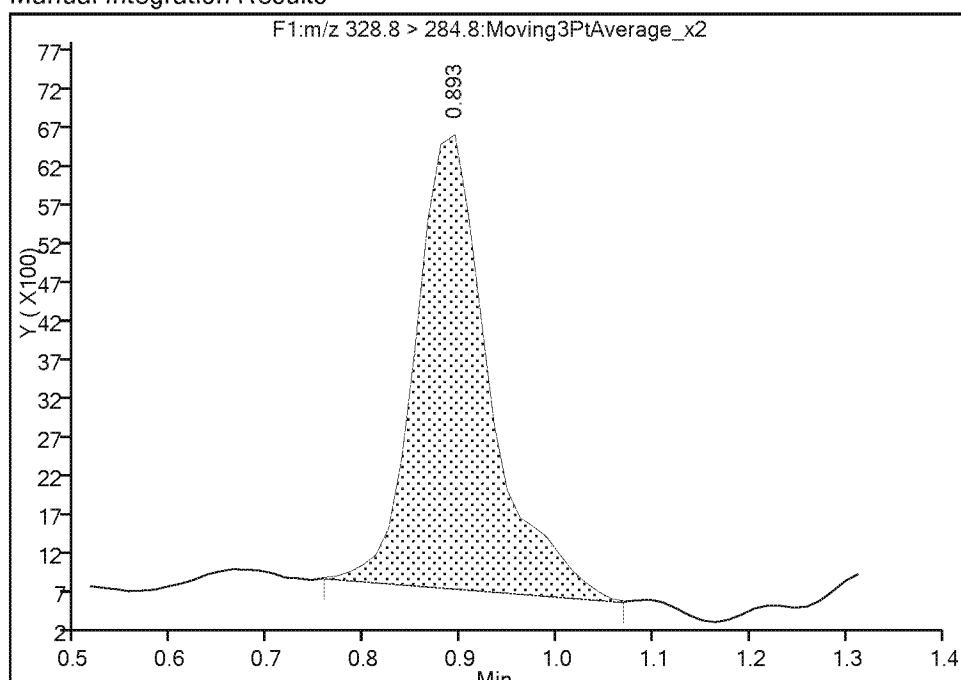

Reviewer: meyera, 10-Oct-2017 11:50:40

Audit Action: Manually Integrated

Audit Reason: Baseline

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10027.d
 Lims ID: std002
 Client ID:
 Sample Type: IC Calib Level: 2
 Inject. Date: 10-Oct-2017 09:38:42 ALS Bottle#: 3 Worklist Smp#: 4
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: L2
 Misc. Info.: HFPO17J10
 Operator ID: JBH Instrument ID: LC_LCMS7
 Sublist: chrom-HFPO*sub1

Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 10-Oct-2017 12:51:46 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK005

First Level Reviewer: meyera Date: 10-Oct-2017 11:50:49

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

* 2 13C3 HFPO-DA (IS)

331.8 > 286.8 0.880 0.880 0.0 745227 10.0 452

\$ 3 13C3 HFPO-DA

331.8 > 286.8 0.880 0.880 0.0 1.000 745227 10.2 452

1 Perfluoro(2-propoxypropanoic) acid

328.8 > 284.8 0.880 0.885 -0.005 1.000 63823 0.6303 36.5

Reagents:

HFPO_CAL-2_00032 Amount Added: 1.00 Units: mL

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20171010-63483.b\\hfp0717J10027.d

Injection Date: 10-Oct-2017 09:38:42 Instrument ID: LC_LCMS7

Lims ID: std002

Client ID:

Operator ID: JBH ALS Bottle#: 3 Worklist Smp#: 4


Injection Vol: 20.0 ul Dil. Factor: 1.0000

Method: HFPO Limit Group: LC - 8321A_HFPO_Du

* 2 13C3 HFPO-DA (IS)

\$ 3 13C3 HFPO-DA

1 Perfluoro(2-propoxypropanoic) acid

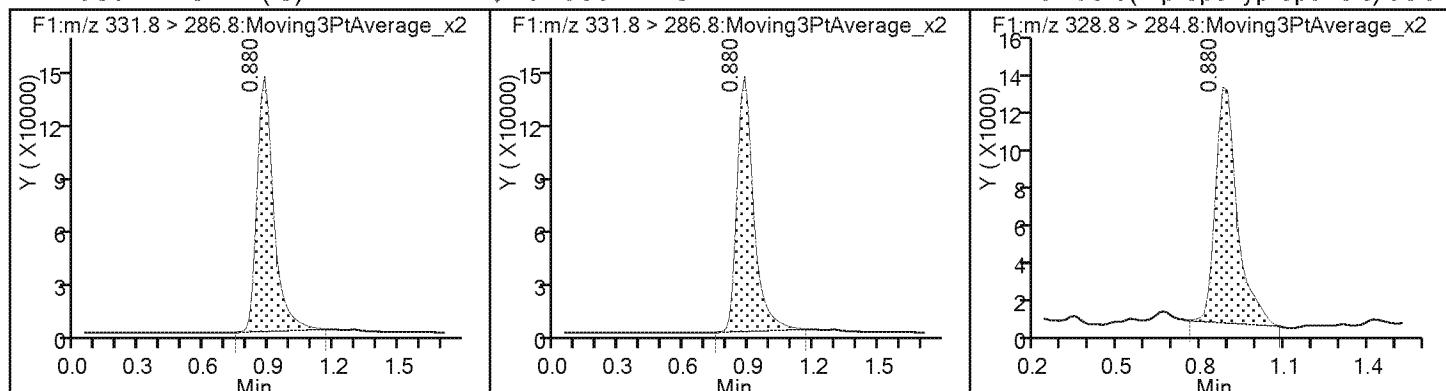

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10028.d
 Lims ID: std003
 Client ID:
 Sample Type: IC Calib Level: 3
 Inject. Date: 10-Oct-2017 09:41:56 ALS Bottle#: 4 Worklist Smp#: 5
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: L3
 Misc. Info.: HFPO17J10
 Operator ID: JBH Instrument ID: LC_LCMS7
 Sublist: chrom-HFPO*sub1

Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 10-Oct-2017 12:51:47 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK005

First Level Reviewer: meyera Date: 10-Oct-2017 11:50:52

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA
 331.8 > 286.8 0.880 0.880 0.0 1.000 750427 10.3 417
 * 2 13C3 HFPO-DA (IS)
 331.8 > 286.8 0.880 0.880 0.0 1.000 750427 10.0 417
 1 Perfluoro(2-propoxypropanoic) acid
 328.8 > 284.8 0.880 0.885 -0.005 1.000 89272 0.9595 50.3

Reagents:

HFPO_CAL-3_00031 Amount Added: 1.00 Units: mL

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20171010-63483.b\\hfp0717J10028.d

Injection Date: 10-Oct-2017 09:41:56

Instrument ID: LC_LCMS7

Lims ID: std003

Client ID:

Operator ID: JBH

ALS Bottle#: 4 Worklist Smp#: 5

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A_HFPO_Du

\$ 3 13C3 HFPO-DA

* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

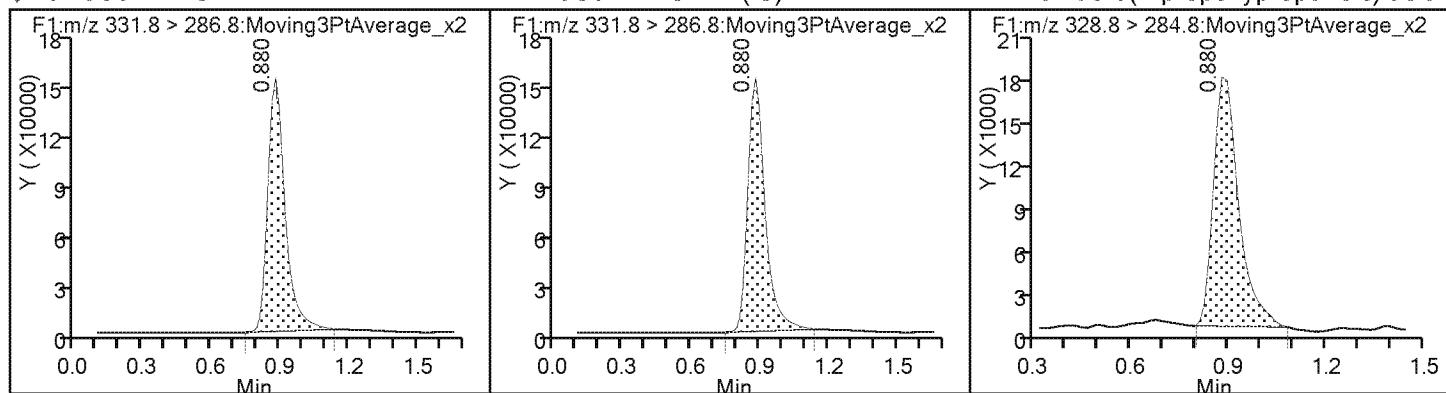

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10029.d
 Lims ID: std004
 Client ID:
 Sample Type: IC Calib Level: 4
 Inject. Date: 10-Oct-2017 09:45:11 ALS Bottle#: 5 Worklist Smp#: 6
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: L4
 Misc. Info.: HFPO17J10
 Operator ID: JBH Instrument ID: LC_LCMS7
 Sublist: chrom-HFPO*sub1

Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 10-Oct-2017 12:51:47 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK005

First Level Reviewer: meyera Date: 10-Oct-2017 11:50:55

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

* 2 13C3 HFPO-DA (IS)

331.8 > 286.8 0.880 0.880 0.0 718028 10.0 438

\$ 3 13C3 HFPO-DA

331.8 > 286.8 0.880 0.880 0.0 1.000 718028 9.82 438

1 Perfluoro(2-propoxypropanoic) acid

328.8 > 284.8 0.880 0.885 -0.005 1.000 167109 2.08 143

Reagents:

HFPO_CAL-4_00031 Amount Added: 1.00 Units: mL

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20171010-63483.b\\hfp0717J10029.d

Injection Date: 10-Oct-2017 09:45:11 Instrument ID: LC_LCMS7

Lims ID: std004

Client ID:

Operator ID: JBH ALS Bottle#: 5 Worklist Smp#: 6


Injection Vol: 20.0 ul Dil. Factor: 1.0000

Method: HFPO Limit Group: LC - 8321A_HFPO_Du

* 2 13C3 HFPO-DA (IS)

\$ 3 13C3 HFPO-DA

1 Perfluoro(2-propoxypropanoic) acid

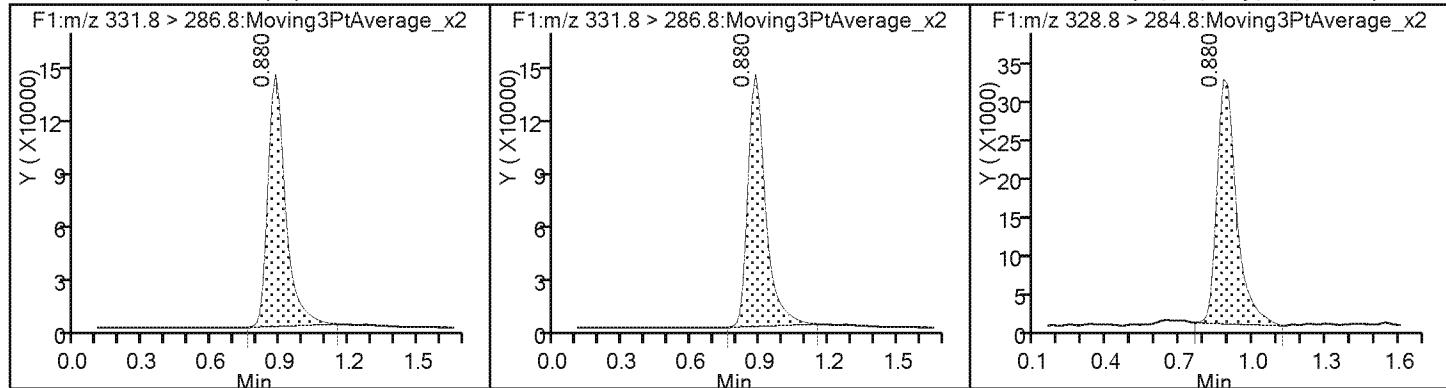

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10030.d
 Lims ID: std005
 Client ID:
 Sample Type: IC Calib Level: 5
 Inject. Date: 10-Oct-2017 09:48:25 ALS Bottle#: 6 Worklist Smp#: 7
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: L5
 Misc. Info.: HFPO17J10
 Operator ID: JBH Instrument ID: LC_LCMS7
 Sublist: chrom-HFPO*sub1

Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 10-Oct-2017 12:51:48 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK005

First Level Reviewer: meyera Date: 10-Oct-2017 11:50:57

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA
331.8 > 286.8 0.880 0.880 0.0 1.000 744600 10.2 433

* 2 13C3 HFPO-DA (IS)
331.8 > 286.8 0.880 0.880 0.0 1.000 744600 10.0 433

1 Perfluoro(2-propoxypropanoic) acid
328.8 > 284.8 0.893 0.885 0.008 1.000 378047 4.80 223

Reagents:

HFPO_CAL-5_00070 Amount Added: 1.00 Units: mL

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20171010-63483.b\\hfpo717J10030.d

Injection Date: 10-Oct-2017 09:48:25 Instrument ID: LC_LCMS7

Lims ID: std005

Client ID:

Operator ID: JBH ALS Bottle#: 6 Worklist Smp#: 7


Injection Vol: 20.0 ul Dil. Factor: 1.0000

Method: HFPO Limit Group: LC - 8321A_HFPO_Du

\$ 3 13C3 HFPO-DA

* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

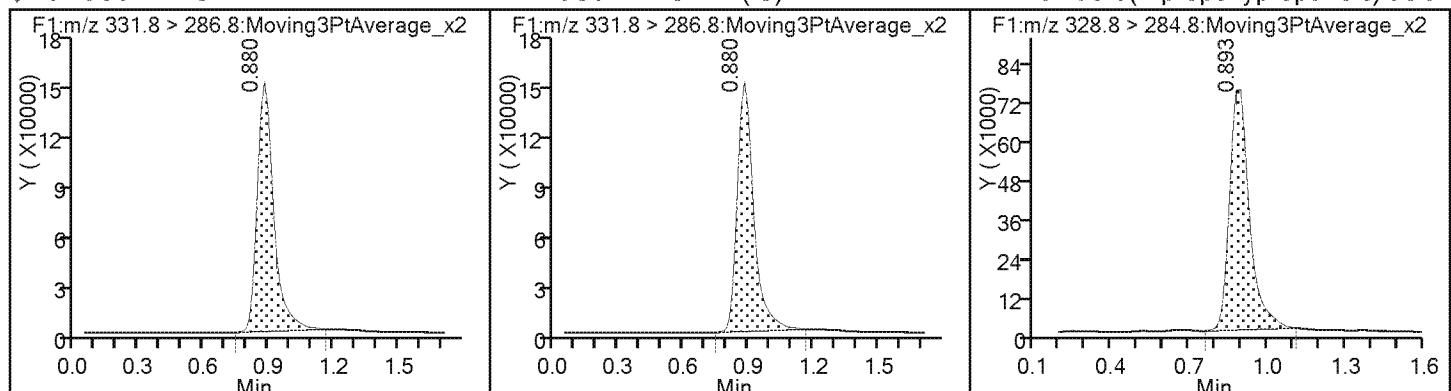

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10031.d
 Lims ID: std006
 Client ID:
 Sample Type: IC Calib Level: 6
 Inject. Date: 10-Oct-2017 09:51:39 ALS Bottle#: 7 Worklist Smp#: 8
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: L6
 Misc. Info.: HFPO17J10
 Operator ID: JBH Instrument ID: LC_LCMS7
 Sublist: chrom-HFPO*sub1

Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 10-Oct-2017 12:51:49 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM
 Process Host: XAWRK005

First Level Reviewer: meyera Date: 10-Oct-2017 11:51:00

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

* 2 13C3 HFPO-DA (IS)
 331.8 > 286.8 0.880 0.880 0.0 731935 10.0 379
 \$ 3 13C3 HFPO-DA
 331.8 > 286.8 0.880 0.880 0.0 1.000 731935 10.0 379
 1 Perfluoro(2-propoxypropanoic) acid
 328.8 > 284.8 0.880 0.885 -0.005 1.000 739399 9.77 298

Reagents:

HFPO_CAL-6_00070 Amount Added: 1.00 Units: mL

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20171010-63483.b\\hfp0717J10031.d

Injection Date: 10-Oct-2017 09:51:39

Instrument ID: LC_LCMS7

Lims ID: std006

Client ID:

Operator ID: JBH

ALS Bottle#: 7 Worklist Smp#: 8

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A_HFPO_Du

* 2 13C3 HFPO-DA (IS)

\$ 3 13C3 HFPO-DA

1 Perfluoro(2-propoxypropanoic) acid

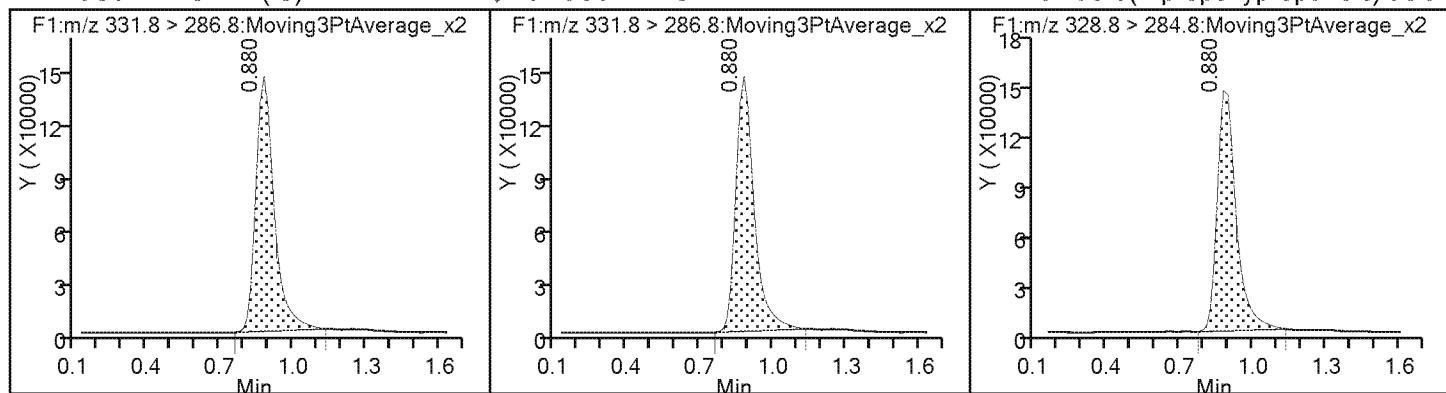

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10032.d
 Lims ID: std007
 Client ID:
 Sample Type: IC Calib Level: 7
 Inject. Date: 10-Oct-2017 09:54:53 ALS Bottle#: 8 Worklist Smp#: 9
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: L7
 Misc. Info.: HFPO17J10
 Operator ID: JBH Instrument ID: LC_LCMS7
 Sublist: chrom-HFPO*sub1

Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 10-Oct-2017 12:51:50 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM
 Process Host: XAWRK005

First Level Reviewer: meyera Date: 10-Oct-2017 11:51:04

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA
 331.8 > 286.8 0.880 0.880 0.0 1.000 729188 9.97 404
 * 2 13C3 HFPO-DA (IS)
 331.8 > 286.8 0.880 0.880 0.0 1.000 729188 10.0 404
 1 Perfluoro(2-propoxypropanoic) acid
 328.8 > 284.8 0.880 0.885 -0.005 1.000 1790812 24.0 386

Reagents:

HFPO_CAL-7_00031 Amount Added: 1.00 Units: mL

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20171010-63483.b\\hfp0717J10032.d

Injection Date: 10-Oct-2017 09:54:53 Instrument ID: LC_LCMS7

Lims ID: std007

Client ID:

Operator ID: JBH

ALS Bottle#: 8 Worklist Smp#: 9

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A_HFPO_Du

\$ 3 13C3 HFPO-DA

* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

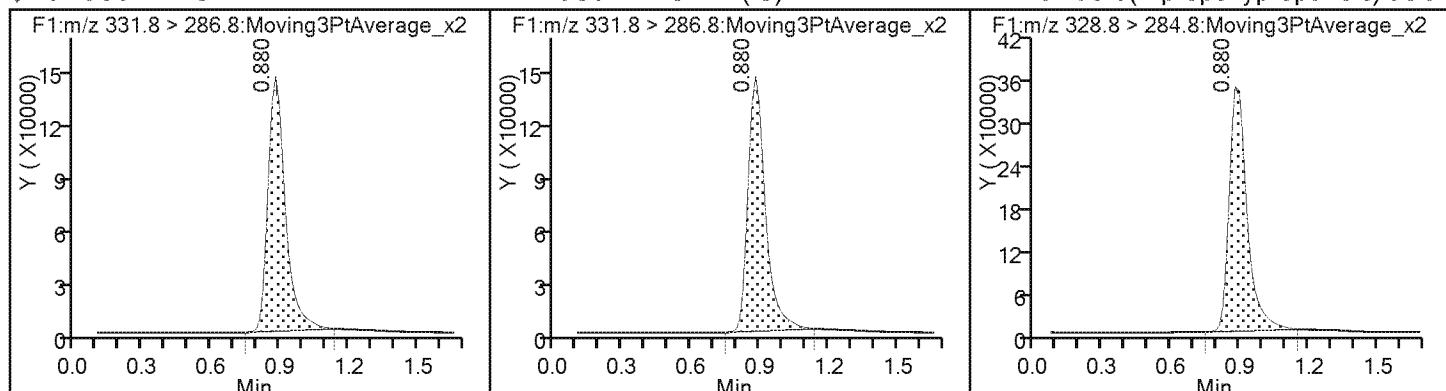

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d
 Lims ID: std008
 Client ID:
 Sample Type: IC Calib Level: 8
 Inject. Date: 10-Oct-2017 09:58:07 ALS Bottle#: 9 Worklist Smp#: 10
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: L8
 Misc. Info.: HFPO17J10
 Operator ID: JBH Instrument ID: LC_LCMS7
 Sublist: chrom-HFPO*sub1

Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 10-Oct-2017 12:51:51 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM
 Process Host: XAWRK005

First Level Reviewer: meyera Date: 10-Oct-2017 11:51:08

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

* 2 13C3 HFPO-DA (IS)
 331.8 > 286.8 0.880 0.880 0.0 701420 10.0 373
 \$ 3 13C3 HFPO-DA
 331.8 > 286.8 0.880 0.880 0.0 1.000 701420 9.59 373
 1 Perfluoro(2-propoxypropanoic) acid
 328.8 > 284.8 0.893 0.885 0.008 1.000 3654104 51.3 421

Reagents:

HFPO_CAL-8_00031 Amount Added: 1.00 Units: mL

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20171010-63483.b\\hfp0717J10033.d

Injection Date: 10-Oct-2017 09:58:07 Instrument ID: LC_LCMS7

Lims ID: std008

Client ID:

Operator ID: JBH

ALS Bottle#: 9 Worklist Smp#: 10

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A_HFPO_Du

* 2 13C3 HFPO-DA (IS)

\$ 3 13C3 HFPO-DA

1 Perfluoro(2-propoxypropanoic) acid

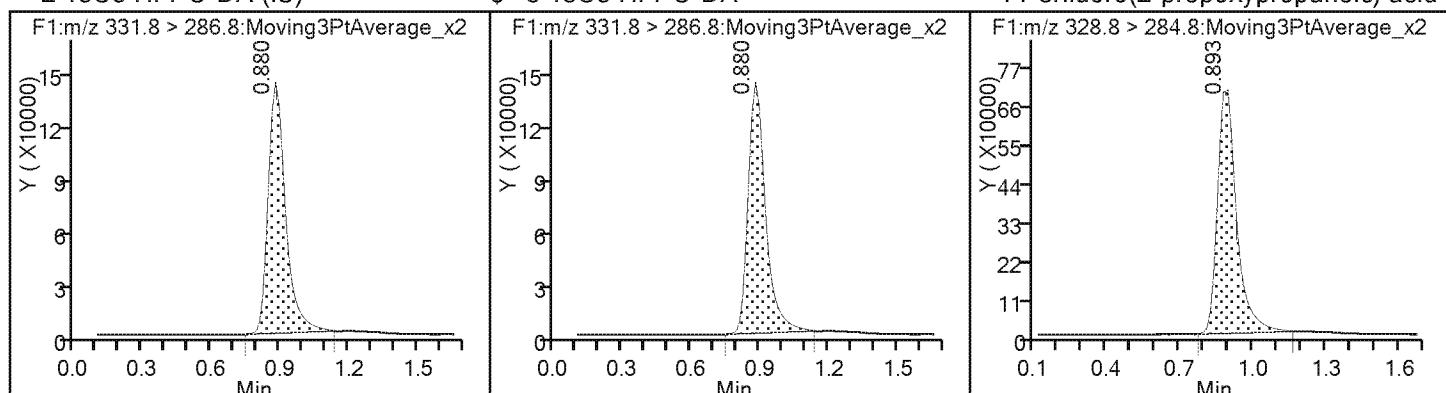

FORM VII
LCMS CONTINUING CALIBRATION DATA

Lab Name: TestAmerica Denver

Job No.: 320-35355-1

SDG No.:

Lab Sample ID: ICV 280-387775/13

Calibration Date: 09/14/2017 15:10

Instrument ID: LC_LCMS7

Calib Start Date: 09/14/2017 14:40

GC Column: Synergi Hydro ID:

Calib End Date: 09/14/2017 15:01

Lab File ID: hfpo717II14062.d

Conc. Units: ug/L

ANALYTE	CURVE TYPE	AVE RRF	RRF	MIN RRF	CALC AMOUNT	SPIKE AMOUNT	%D	MAX %D
Perfluoro(2-propoxypropanoic acid	Lin1		0.9462		1.89	2.00	-5.3	20.0
13C3 HFPO-DA	Ave	192740	197806		10.3	10.0	2.6	

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20170915-62647.b\hfpo717I14062.d
 Lims ID: ICV
 Client ID:
 Sample Type: ICV
 Inject. Date: 14-Sep-2017 15:10:31 ALS Bottle#: 10 Worklist Smp#: 13
 Injection Vol: 10.0 ul Dil. Factor: 1.0000
 Sample Info: ICV
 Misc. Info.: HFPO17I14
 Operator ID: JBH Instrument ID: LC_LCMS7
 Sublist:
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20170915-62647.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 15-Sep-2017 07:29:44 Calib Date: 14-Sep-2017 15:01:22
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20170915-62647.b\hfpo717I14059.d
 Column 1 : Det: F1:MRM
 Process Host: XAWRK034

First Level Reviewer: meyera Date: 15-Sep-2017 07:28:43

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA
 331.8 > 286.8 0.988 0.981 0.007 1.000 1978058 10.3 436
 * 2 13C3 HFPO-DA (IS)
 331.8 > 286.8 0.988 0.981 0.007 1.000 1978058 10.0 436
 1 Perfluoro(2-propoxypropanoic) acid
 328.8 > 284.8 0.988 0.986 0.002 1.000 374307 1.89 162

Reagents:

HFPO_ICV_00031 Amount Added: 1.00 Units: mL

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20170915-62647.b\\hfp0717\\14062.d

Injection Date: 14-Sep-2017 15:10:31 Instrument ID: LC_LCMS7

Lims ID: ICV

Client ID:

Operator ID: JBH ALS Bottle#: 10 Worklist Smp#: 13


Injection Vol: 10.0 ul Dil. Factor: 1.0000

Method: HFPO Limit Group: LC - 8321A_HFPO_Du

\$ 3 13C3 HFPO-DA

* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

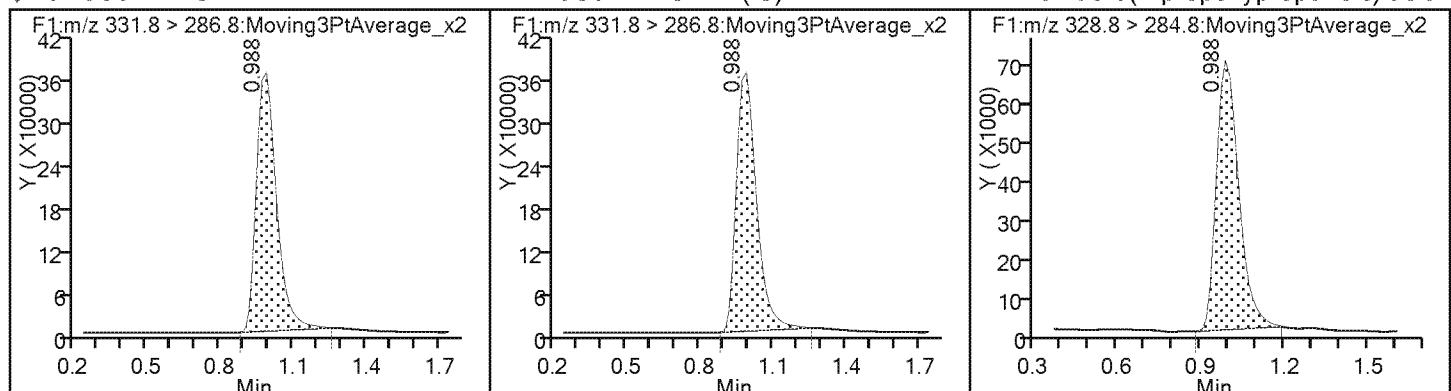

FORM VII
LCMS CONTINUING CALIBRATION DATA

Lab Name: TestAmerica Denver

Job No.: 320-35355-1

SDG No.:

Lab Sample ID: ICV 280-390728/13

Calibration Date: 10/10/2017 10:07

Instrument ID: LC_LCMS7

Calib Start Date: 10/10/2017 09:35

GC Column: Synergi Hydro ID:

Calib End Date: 10/10/2017 09:58

Lab File ID: hfpo717J10036.d

Conc. Units: ug/L

ANALYTE	CURVE TYPE	AVE RRF	RRF	MIN RRF	CALC AMOUNT	SPIKE AMOUNT	%D	MAX %D
HFPO-DA	Lin1		1.154		2.07	2.00	3.3	20.0
13C3 HFPO-DA	Ave	73145	72923		9.97	10.0	-0.3	

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10036.d
 Lims ID: ICV
 Client ID:
 Sample Type: ICV
 Inject. Date: 10-Oct-2017 10:07:48 ALS Bottle#: 10 Worklist Smp#: 13
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: ICV
 Misc. Info.: HFPO17J10
 Operator ID: JBH Instrument ID: LC_LCMS7
 Sublist:
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 10-Oct-2017 12:51:53 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK005

First Level Reviewer: meyera Date: 10-Oct-2017 11:51:34

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA
 331.8 > 286.8 0.880 0.880 0.0 1.000 729225 9.97 396
 * 2 13C3 HFPO-DA (IS)
 331.8 > 286.8 0.880 0.880 0.0 1.000 729225 10.0 396
 1 Perfluoro(2-propoxypropanoic) acid
 328.8 > 284.8 0.893 0.885 0.008 1.000 168368 2.07 111

Reagents:

HFPO_ICV_00032 Amount Added: 1.00 Units: mL

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20171010-63483.b\\hfp0717J10036.d

Injection Date: 10-Oct-2017 10:07:48 Instrument ID: LC_LCMS7

Lims ID: ICV

Client ID:

Operator ID: JBH

ALS Bottle#: 10 Worklist Smp#: 13

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A_HFPO_Du

\$ 3 13C3 HFPO-DA

* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

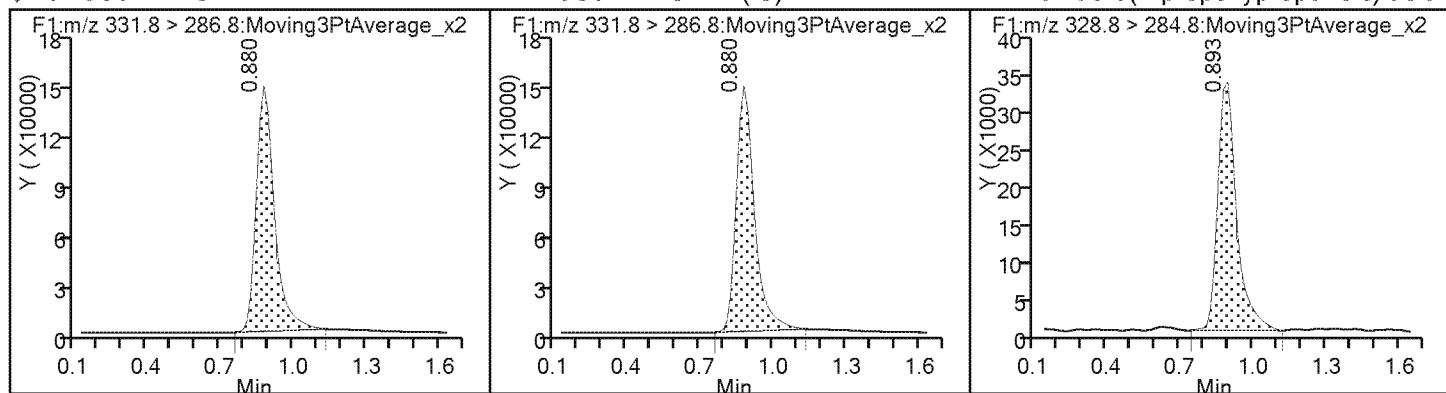

FORM VII
LCMS CONTINUING CALIBRATION DATA

Lab Name: TestAmerica Denver

Job No.: 320-35355-1

SDG No.:

Lab Sample ID: CCV 280-390728/24

Calibration Date: 10/10/2017 10:43

Instrument ID: LC_LCMS7

Calib Start Date: 10/10/2017 09:35

GC Column: Synergi Hydro ID:

Calib End Date: 10/10/2017 09:58

Lab File ID: hfpo717J10047.d

Conc. Units: ug/L

ANALYTE	CURVE TYPE	AVE RRF	RRF	MIN RRF	CALC AMOUNT	SPIKE AMOUNT	%D	MAX %D
HFPO-DA	Lin1		1.012		9.78	10.0	-2.2	20.0
13C3 HFPO-DA	Ave	73145	68787		9.40	10.0	-6.0	

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10047.d
 Lims ID: CCV L6
 Client ID:
 Sample Type: CCV
 Inject. Date: 10-Oct-2017 10:43:29 ALS Bottle#: 7 Worklist Smp#: 24
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: CCV L6
 Misc. Info.: HFPO17J10
 Operator ID: JBH Instrument ID: LC_LCMS7
 Sublist: chrom-HFPO*sub1
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 10-Oct-2017 12:52:02 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICAL File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d
 Column 1 : Det: F1:MRM
 Process Host: XAWRK005

First Level Reviewer: meyera Date: 10-Oct-2017 11:52:05

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA
 331.8 > 286.8 0.839 0.880 -0.041 1.000 687867 9.40 327
 * 2 13C3 HFPO-DA (IS)
 331.8 > 286.8 0.839 0.880 -0.041 687867 10.0 327
 1 Perfluoro(2-propoxypropanoic) acid
 328.8 > 284.8 0.839 0.885 -0.046 1.000 696191 9.78 224

Reagents:

HFPO_CAL-6_00070 Amount Added: 1.00 Units: mL

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20171010-63483.b\\hfp0717J10047.d

Injection Date: 10-Oct-2017 10:43:29 Instrument ID: LC_LCMS7

Lims ID: CCV L6

Client ID:

Operator ID: JBH

ALS Bottle#: 7 Worklist Smp#: 24

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A_HFPO_Du

\$ 3 13C3 HFPO-DA

* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

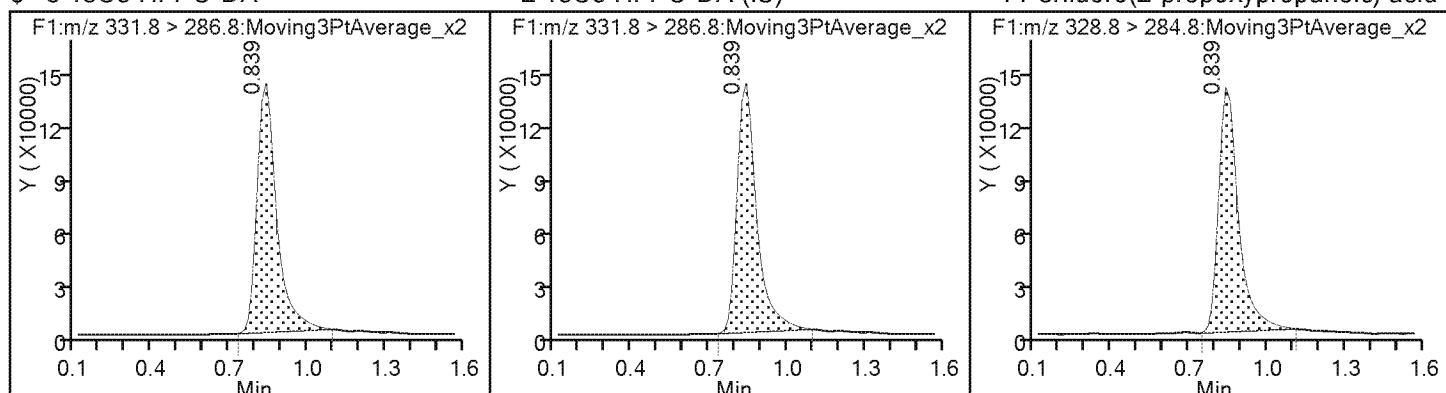

FORM VII
LCMS CONTINUING CALIBRATION DATA

Lab Name: TestAmerica Denver

Job No.: 320-35355-1

SDG No.:

Lab Sample ID: CCV 280-403246/3

Calibration Date: 01/30/2018 09:12

Instrument ID: LC_LCMS7

Calib Start Date: 10/10/2017 09:35

GC Column: Synergi Hydro ID:

Calib End Date: 10/10/2017 09:58

Lab File ID: hfpo718A30003.d

Conc. Units: ug/L

ANALYTE	CURVE TYPE	AVE RRF	RRF	MIN RRF	CALC AMOUNT	SPIKE AMOUNT	%D	MAX %D
HFPO-DA	Lin1		0.9303		8.98	10.0	-10.2	20.0
13C3 HFPO-DA	Ave	73145	76648		10.5	10.0	4.8	

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30003.d
 Lims ID: CCV L6
 Client ID:
 Sample Type: CCV
 Inject. Date: 30-Jan-2018 09:12:33 ALS Bottle#: 7 Worklist Smp#: 3
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: CCV L6
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Sublist: chrom-HFPO*sub1

Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 30-Jan-2018 14:49:56 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICAL File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM
 Process Host: XAWRK010

First Level Reviewer: meyera Date: 30-Jan-2018 14:47:10

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA
 331.8 > 286.8 1.042 0.961 0.081 1.000 766482 10.5 1557
 * 2 13C3 HFPO-DA (IS)
 331.8 > 286.8 1.042 0.961 0.081 1.000 766482 10.0 1557
 1 Perfluoro(2-propoxypropanoic) acid
 328.8 > 284.8 1.042 0.988 0.054 1.000 713028 8.98 144

Reagents:

HFPO_CAL-6_00079 Amount Added: 1.00 Units: mL

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20180130-66815.b\\hfpo718A30003.d

Injection Date: 30-Jan-2018 09:12:33

Instrument ID: LC_LCMS7

Lims ID: CCV L6

Client ID:

Operator ID: JBH

ALS Bottle#: 7 Worklist Smp#: 3

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A_HFPO_Du

\$ 3 13C3 HFPO-DA

* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

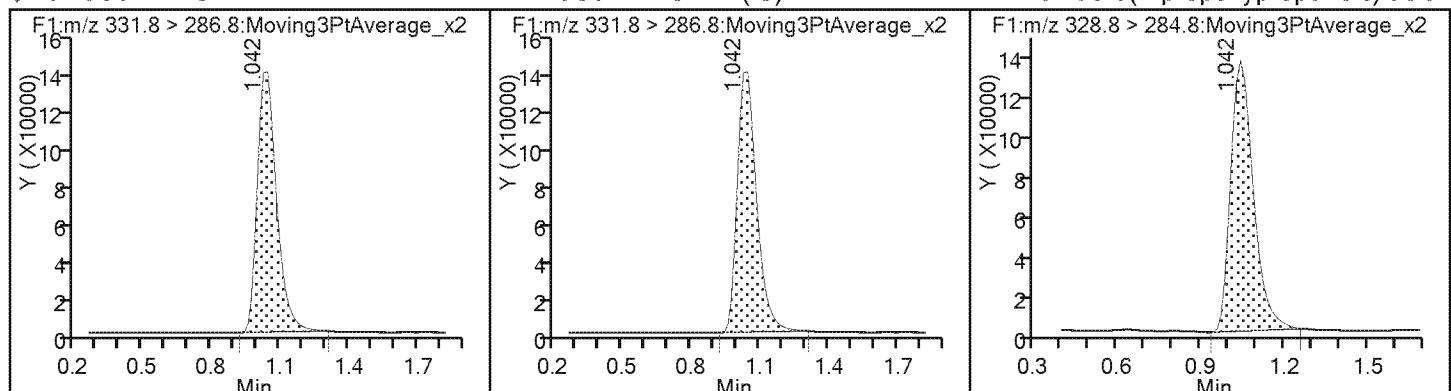

FORM VII
LCMS CONTINUING CALIBRATION DATA

Lab Name: TestAmerica Denver Job No.: 320-35355-1
SDG No.:
Lab Sample ID: CCV 280-403246/12 Calibration Date: 01/30/2018 09:41
Instrument ID: LC_LCMS7 Calib Start Date: 10/10/2017 09:35
GC Column: Synergi Hydro ID: Calib End Date: 10/10/2017 09:58
Lab File ID: hfpo718A30012.d Conc. Units: ug/L

ANALYTE	CURVE TYPE	AVE RRF	RRF	MIN RRF	CALC AMOUNT	SPIKE AMOUNT	%D	MAX %D
HFPO-DA	Lin1		0.9119		4.29	5.00	-14.2	20.0
13C3 HFPO-DA	Ave	73145	102446		14.0	10.0	40.1	

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30012.d
 Lims ID: CCV L5
 Client ID:
 Sample Type: CCV
 Inject. Date: 30-Jan-2018 09:41:47 ALS Bottle#: 6 Worklist Smp#: 12
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: CCV L5
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Sublist: chrom-HFPO*sub1

Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 30-Jan-2018 14:50:04 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK010

First Level Reviewer: meyera Date: 30-Jan-2018 14:48:02

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 > 286.8 0.961 0.961 0.0 1.000 1024461 14.0 1765

* 2 13C3 HFPO-DA (IS)

331.8 > 286.8 0.961 0.961 0.0 1.000 1024461 10.0 1765

1 Perfluoro(2-propoxypropanoic) acid

328.8 > 284.8 0.961 0.988 -0.027 1.000 467090 4.29 226

Reagents:

HFPO_CAL-5_00079 Amount Added: 1.00 Units: mL

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20180130-66815.b\\hfpo718A30012.d

Injection Date: 30-Jan-2018 09:41:47

Instrument ID: LC_LCMS7

Lims ID: CCV L5

Client ID:

Operator ID: JBH

ALS Bottle#: 6 Worklist Smp#: 12

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A_HFPO_Du

\$ 3 13C3 HFPO-DA

* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

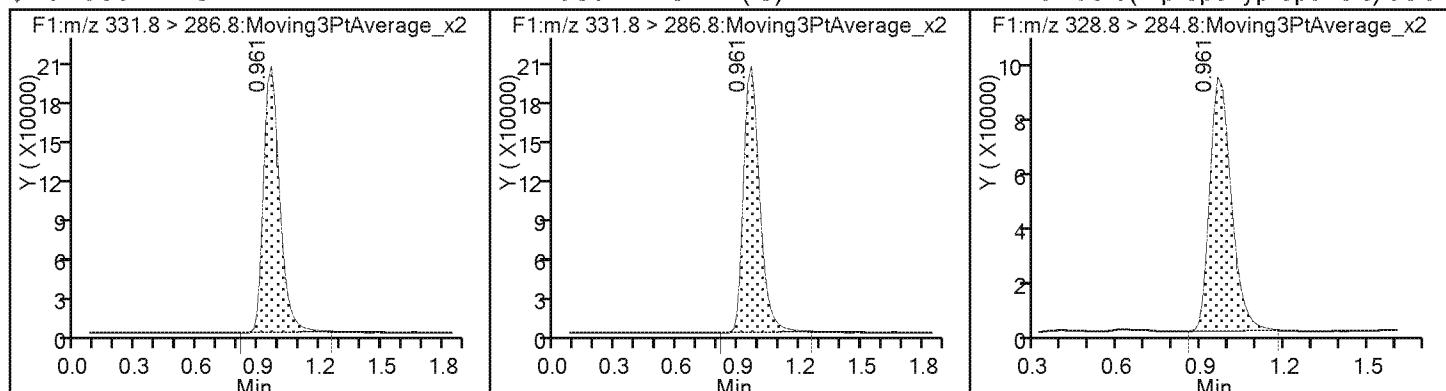

FORM VII
LCMS CONTINUING CALIBRATION DATA

Lab Name: TestAmerica Denver Job No.: 320-35355-1
SDG No.:
Lab Sample ID: CCV 280-403246/21 Calibration Date: 01/30/2018 10:11
Instrument ID: LC_LCMS7 Calib Start Date: 10/10/2017 09:35
GC Column: Synergi Hydro ID: Calib End Date: 10/10/2017 09:58
Lab File ID: hfpo718A30021.d Conc. Units: ug/L

ANALYTE	CURVE TYPE	AVE RRF	RRF	MIN RRF	CALC AMOUNT	SPIKE AMOUNT	%D	MAX %D
HFPO-DA	Lin1		0.9411		9.08	10.0	-9.2	20.0
13C3 HFPO-DA	Ave	73145	98761		13.5	10.0	35.0	

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30021.d
 Lims ID: CCV L6
 Client ID:
 Sample Type: CCV
 Inject. Date: 30-Jan-2018 10:11:04 ALS Bottle#: 7 Worklist Smp#: 21
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: CCV L6
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Sublist: chrom-HFPO*sub1
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 30-Jan-2018 14:50:10 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK010

First Level Reviewer: meyera Date: 30-Jan-2018 14:49:13

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 > 286.8 0.961 0.961 0.0 1.000 987612 13.5 1456

* 2 13C3 HFPO-DA (IS)

331.8 > 286.8 0.961 0.961 0.0 1.000 987612 10.0 1456

1 Perfluoro(2-propoxypropanoic) acid

328.8 > 284.8 0.961 0.988 -0.027 1.000 929479 9.08 396

Reagents:

HFPO_CAL-6_00079 Amount Added: 1.00 Units: mL

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20180130-66815.b\\hfp0718A30021.d

Injection Date: 30-Jan-2018 10:11:04

Instrument ID: LC_LCMS7

Lims ID: CCV L6

Client ID:

Operator ID: JBH

ALS Bottle#: 7 Worklist Smp#: 21

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A_HFPO_Du

\$ 3 13C3 HFPO-DA

* 2 13C3 HFPO-DA (IS)


1 Perfluoro(2-propoxypropanoic) acid

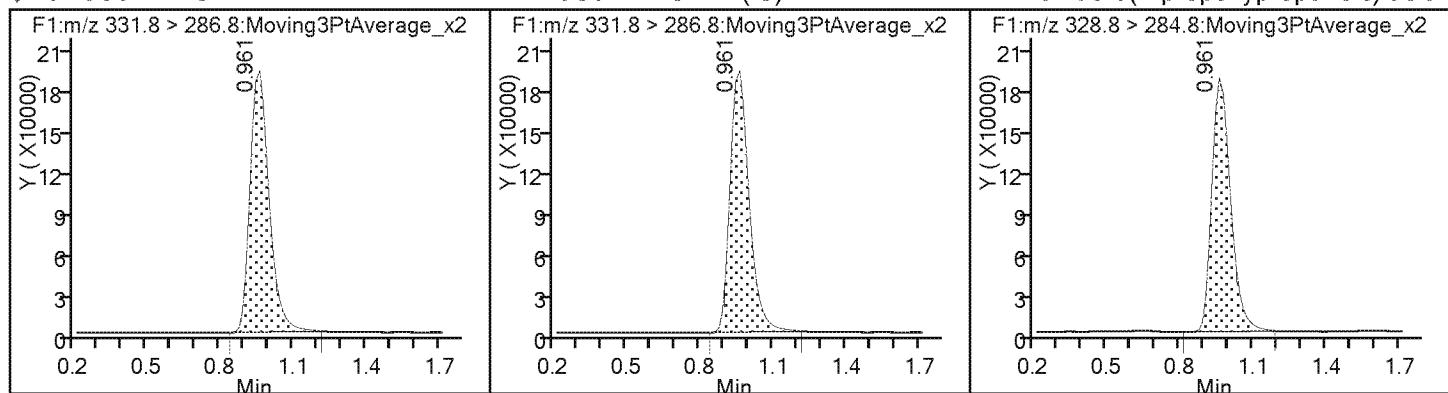

File: C:\MassLynx8321.PRO\ACQ\QDB\HFPO\MRM.lpr

Instrument: XEVO-TQMS\PVBA453

Printed: Tuesday, January 30, 2018 08:19:42 Mountain Standard Time

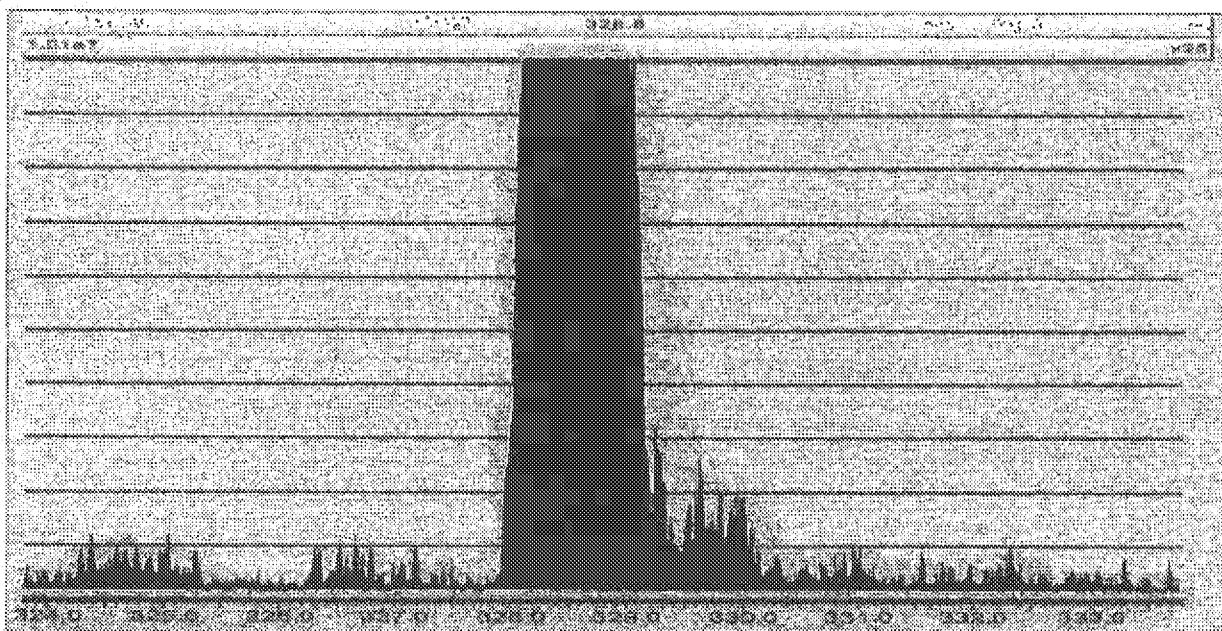

Type	Start Mass	End Mass	Set Mass
MG1 Scan	323.80	333.80	
Sources (ES-)	Settings	Readbacks	
Capillary (kV)	0.50	0.54	
Cone (V)	10.00	-21.06	
Extractor (V)	3.00	-10.61	
Source Temperature (°C)	120	120	
Desolvation Temperature (°C)	200	200	
Cone Gas Flow (L/Hr)	50	50	
Desolvation Gas Flow (L/Hr)	800	791	
Collision Gas Flow (mL/Min)	0.15	0.04	
Analyser	Settings	Readbacks	
LM 1 Resolution	2.8		
HM 1 Resolution	14.8		
Ion Energy 1	0.7		
MS Mode Collision Energy	7.00		
MSMS Mode Collision Energy	20.00		
MS Mode Entrance	0.60		
MS Mode Exit	0.60		
Gas On MS Mode Entrance	0.60		
Gas On MS Mode Exit	0.60		
Gas On MSMS Mode Entrance	0.60		
Gas On MSMS Mode Exit	0.60		
Gas Off MS Mode Entrance	30.00		
Gas Off MS Mode Exit	30.00		
Gas Off MSMS Mode Entrance	2.00		
Gas Off MSMS Mode Exit	2.00		
ScanWave MS Mode Entrance	0.60		
ScanWave MS Mode Exit	0.60		
ScanWave MSMS Mode Entrance	0.60		
ScanWave MSMS Mode Exit	0.60		
LM 2 Resolution	2.9		
HM 2 Resolution	14.7		
Ion Energy 2	0.3		

Thangasim P. 02/02/18

File: C:\MassLynx\8321.PROVACQUDBVHFPOMRM.ipr

Instrument: XEVO-TQMSIVBA463

Printed: Tuesday, January 30, 2018 08:19:42 Mountain Standard Time

Multiplier 623.57
Active Reservoir A

Pressure Gauges
Collision Cell Pressure (mbar) 7.878782e-006

Instrument Configuration**Automatic Mode**

MS Inter-scan delay (secs) 0.005
Polarity/Mode switch Inter-scan delay (secs) 0.020
Enhanced Inter-scan delay (secs) 0.020
Inter-channel delay - See Tables

MS 1 Delay Table:

R	delay
<= 0.600	0.005
<= 2.000	0.008
<= 4.000	0.010
<= 11.000	0.012
> 11.000	0.014

MS 2 Delay Table:


R	delay
<= 8.000	0.005
<= 25.000	0.005
> 25.000	0.007

Thompson, R 03/02/18

File: C:\MassLynx\8321,PROVACQUDBVFFPOMRM.lpr

Instrument: XEVO-TQMS\PBA463

Printed: Tuesday, January 30, 2018 08:20:26 Mountain Standard Time

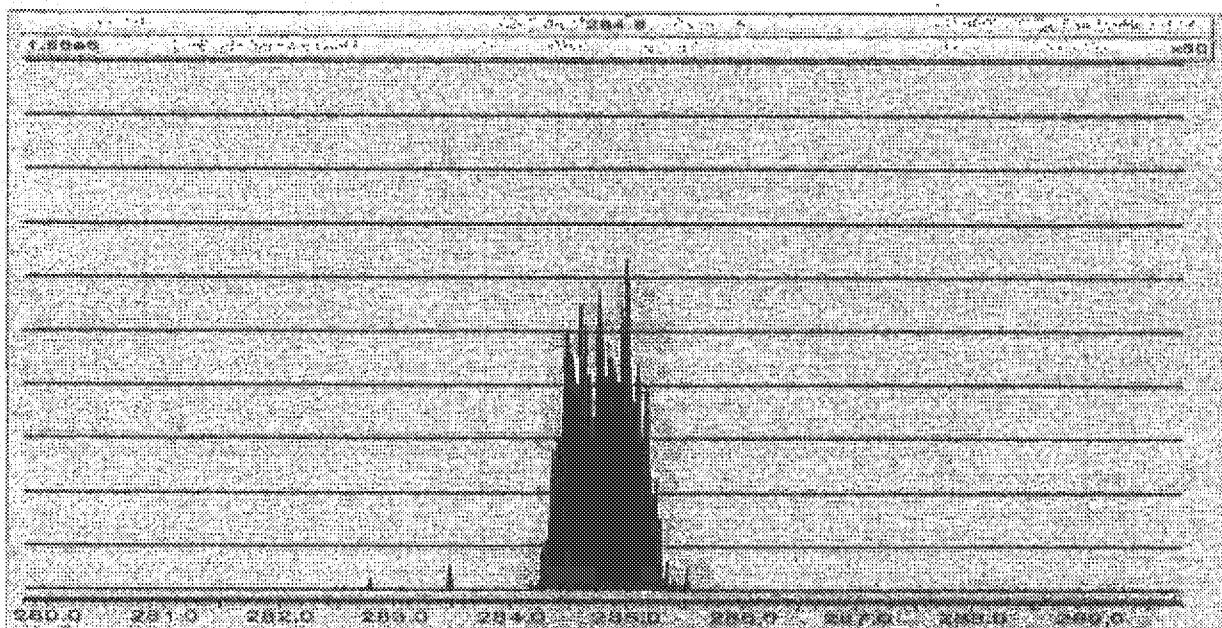

Type	Start Mass	End Mass	Set Mass
Daughter Scan	278.80	289.80	328.80

Source (ES-)	Settings	Readbacks
Capillary (kV)	0.50	0.52
Cone (V)	10.00	-21.06
Extractor (V)	3.00	-10.61
Source Temperature (°C)	120	120
Desolvation Temperature (°C)	200	200
Cone Gas Flow (L/Hr)	50	51
Desolvation Gas Flow (L/Hr)	800	793
Collision Gas Flow (mL/Min)	0.15	0.14

Analyser	Settings	Readbacks
LM 1 Resolution	2.8	
HM 1 Resolution	14.8	
Ion Energy 1	0.7	
MS Moda Collision Energy	7.00	
MSMS Moda Collision Energy	20.00	
MS Moda Entrance	0.60	
MS Moda Exit	0.60	
Gas On MS Moda Entrance	0.60	
Gas On MS Moda Exit	0.60	
Gas On MSMS Moda Entrance	0.60	
Gas On MSMS Moda Exit	0.60	
Gas Off MS Moda Entrance	30.00	
Gas Off MS Moda Exit	30.00	
Gas Off MSMS Moda Entrance	2.00	
Gas Off MSMS Moda Exit	2.00	
ScanWave MS Moda Entrance	0.60	
ScanWave MS Moda Exit	0.60	
ScanWave MSMS Moda Entrance	0.60	
ScanWave MSMS Moda Exit	0.60	
LM 2 Resolution	2.8	
HM 2 Resolution	14.7	
Ion Energy 2	0.3	

Thompson R. 02/02/18

File: C:\MassLynx\8321.PRO\VACQUDBHFPMRM.lpr

Instrument: XEVO-TQMS\IVBA463

Printed: Tuesday, January 30, 2018 08:20:26 Mountain Standard Time

Multiplier 523.81
Active Reservoir A

Pressure Gauges
Collision Cell Pressure (mbar) 1.206243e-003

Instrument Configuration**Automatic Mode**

MS Inter-scan delay (secs) 0.005
Polarity/Mode switch Inter-scan delay (secs) 0.020
Enhanced Inter-scan delay (secs) 0.020
Inter-channel delay - See Tables

MS 1 Delay Table:

R	delay
<= 0.600	0.006
<= 2.000	0.008
<= 4.000	0.010
<= 11.000	0.012
> 11.000	0.014

MS 2 Delay Table:

R	delay
<= 8.000	0.005
<= 25.000	0.005
> 25.000	0.007

Thomassen P. 02/21/08

File: c:\masslynx\8321.protocqudb\hfpo.exp

Printed: Tuesday, January 30, 2018 08:22:19 Mountain Standard Time

Creation Time Fri 18 Nov 2016 09:08:40
Instrument Identifier XEVO-TQMS/VP/BA453
Version Number 1.0
Duration (min) 2.0
Calibration Filename C:\MassLynx\Intellistart\Results\Unit Mass Resolution\Calibration_20100811
2.cal
Solvent Delay Divert Valve Enabled 0
Number Of Functions 1

Function 1 : MRM of 2 mass pairs, Time 0.00 to 2.00, ES-

Type MRM
Ion Mode ES-
Inter Channel Delay (sec) -1.000
InterScan Time (sec) -1.000
Span (Da) 0.6
Start Time (min) 0.0
End Time (min) 2.0
Ch Pcnt(Da) Dau(Da) Dwell(s) Cone(V) Coll(eV) Delay(s) Compound
1 328.80 284.80 0.400 10.00 7.00 -1.000 HFPO
2 331.80 286.80 0.400 10.00 7.00 -1.000 HFPO IS

Thronecell P. cube/p

FORM I
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: TestAmerica Denver Job No.: 320-35355-1

SDG No.: _____

Client Sample ID: _____ Lab Sample ID: MB 280-403139/1-A

Matrix: Water Lab File ID: hfpo718A30004.d

Analysis Method: 8321A Date Collected: _____

Extraction Method: 3535 Date Extracted: 01/29/2018 21:37

Sample wt/vol: 250 (mL) Date Analyzed: 01/30/2018 09:15

Con. Extract Vol.: 5 (mL) Dilution Factor: 1

Injection Volume: 20 (uL) GC Column: Synergi Hydro ID: _____

% Moisture: GPC Cleanup: (Y/N) N

Analysis Batch No.: 403246 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	<0.010		0.010	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	111		50-200

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30004.d
 Lims ID: MB 280-403139/1-A
 Client ID:
 Sample Type: MB
 Inject. Date: 30-Jan-2018 09:15:47 ALS Bottle#: 11 Worklist Smp#: 4
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: MB280-403139/1-A
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 30-Jan-2018 14:49:56 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK010

First Level Reviewer: meyera Date: 30-Jan-2018 14:47:14

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA
 331.8 > 286.8 0.988 0.961 0.027 1.000 808587 11.1 2714
 * 2 13C3 HFPO-DA (IS)
 331.8 > 286.8 0.988 0.961 0.027 808587 10.0 2714

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20180130-66815.b\\hfpo718A30004.d

Injection Date: 30-Jan-2018 09:15:47

Instrument ID: LC_LCMS7

Lims ID: MB 280-403139/1-A

Client ID:

Operator ID: JBH

ALS Bottle#: 11 Worklist Smp#: 4

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A_HFPO_Du

\$ 3 13C3 HFPO-DA

* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid (ND)

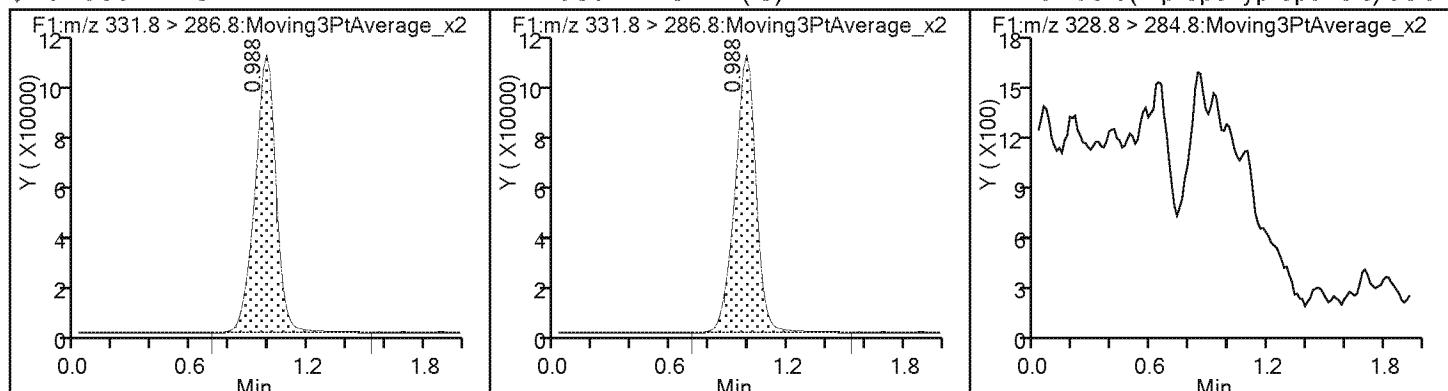

TestAmerica Denver
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30004.d
 Lims ID: MB 280-403139/1-A
 Client ID:
 Sample Type: MB
 Inject. Date: 30-Jan-2018 09:15:47 ALS Bottle#: 11 Worklist Smp#: 4
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: MB280-403139/1-A
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 30-Jan-2018 14:49:56 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK010

First Level Reviewer: meyera Date: 30-Jan-2018 14:47:14

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	11.1	110.55

FORM I
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: TestAmerica Denver Job No.: 320-35355-1
SDG No.:
Client Sample ID: Lab Sample ID: ICB 280-390728/11
Matrix: Water Lab File ID: hfpo717J10034.d
Analysis Method: 8321A Date Collected:
Extraction Method: Date Extracted:
Sample wt/vol: 1 (mL) Date Analyzed: 10/10/2017 10:01
Con. Extract Vol.: Dilution Factor: 1
Injection Volume: 20 (uL) GC Column: Synergi Hydro ID:
% Moisture: GPC Cleanup: (Y/N) N
Analysis Batch No.: 390728 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	<0.50		0.50	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	100		50-200

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10034.d
 Lims ID: ICB
 Client ID:
 Sample Type: ICB
 Inject. Date: 10-Oct-2017 10:01:21 ALS Bottle#: 1 Worklist Smp#: 11
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: ICB
 Misc. Info.: HFPO17J10
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 10-Oct-2017 12:51:51 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK005

First Level Reviewer: meyera Date: 10-Oct-2017 11:51:12

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA
 331.8 > 286.8 0.880 0.880 0.0 1.000 732194 10.0 425
 * 2 13C3 HFPO-DA (IS)
 331.8 > 286.8 0.880 0.880 0.0 1.000 732194 10.0 425
 1 Perfluoro(2-propoxypropanoic) acid
 328.8 > 284.8 0.880 0.885 -0.005 1.000 13993 -0.0270 8.1

Reagents:

HFPO_CAL-0_00031 Amount Added: 1.00 Units: mL

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20171010-63483.b\\hfp0717J10034.d

Injection Date: 10-Oct-2017 10:01:21 Instrument ID: LC_LCMS7

Lims ID: ICB

Client ID:

Operator ID: JBH

ALS Bottle#: 1 Worklist Smp#: 11

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A_HFPO_Du

\$ 3 13C3 HFPO-DA

* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

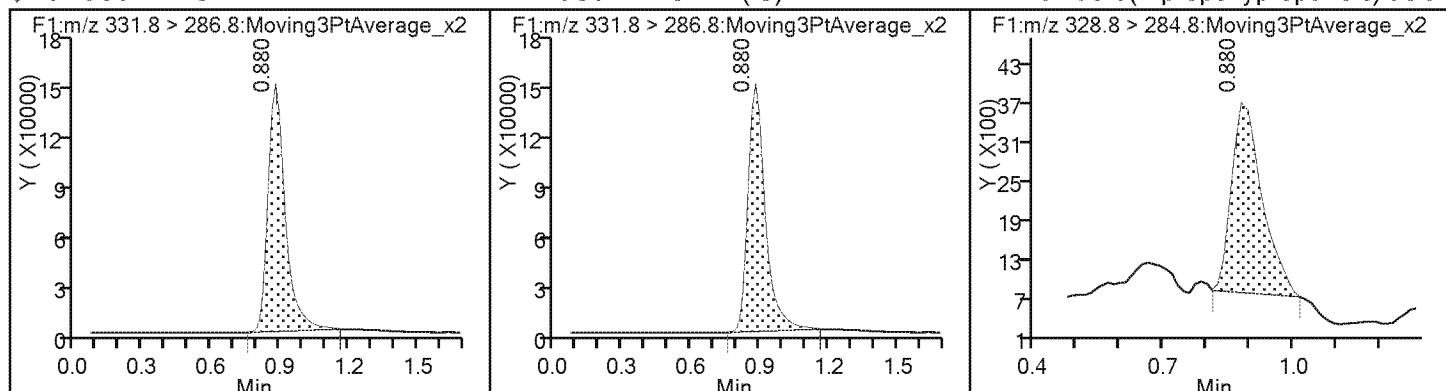

TestAmerica Denver
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10034.d
 Lims ID: ICB
 Client ID:
 Sample Type: ICB
 Inject. Date: 10-Oct-2017 10:01:21 ALS Bottle#: 1 Worklist Smp#: 11
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: ICB
 Misc. Info.: HFPO17J10
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 10-Oct-2017 12:51:51 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK005

First Level Reviewer: meyera Date: 10-Oct-2017 11:51:12

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	10.0	100.10

FORM I
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: TestAmerica Denver Job No.: 320-35355-1

SDG No.: _____

Client Sample ID: _____ Lab Sample ID: LCS 280-403139/2-A

Matrix: Water Lab File ID: hfpo718A30005.d

Analysis Method: 8321A Date Collected: _____

Extraction Method: 3535 Date Extracted: 01/29/2018 21:37

Sample wt/vol: 250 (mL) Date Analyzed: 01/30/2018 09:19

Con. Extract Vol.: 5 (mL) Dilution Factor: 1

Injection Volume: 20 (uL) GC Column: Synergi Hydro ID: _____

% Moisture: GPC Cleanup: (Y/N) N

Analysis Batch No.: 403246 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	0.225		0.010	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	107		50-200

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30005.d
 Lims ID: LCS 280-403139/2-A
 Client ID:
 Sample Type: LCS
 Inject. Date: 30-Jan-2018 09:19:02 ALS Bottle#: 12 Worklist Smp#: 5
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: LCS280-403139/2-A
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 30-Jan-2018 14:49:56 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK010

First Level Reviewer: meyera Date: 30-Jan-2018 14:47:19

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA
 331.8 > 286.8 0.934 0.961 -0.027 1.000 781298 10.7 1973
 * 2 13C3 HFPO-DA (IS)
 331.8 > 286.8 0.934 0.961 -0.027 1.000 781298 10.0 1973
 1 Perfluoro(2-propoxypropanoic) acid
 328.8 > 284.8 0.934 0.988 -0.054 1.000 906720 11.3 457

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20180130-66815.b\\hfpo718A30005.d

Injection Date: 30-Jan-2018 09:19:02

Instrument ID: LC_LCMS7

Lims ID: LCS 280-403139/2-A

Client ID:

Operator ID: JBH

ALS Bottle#: 12 Worklist Smp#: 5

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A_HFPO_Du

\$ 3 13C3 HFPO-DA

* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

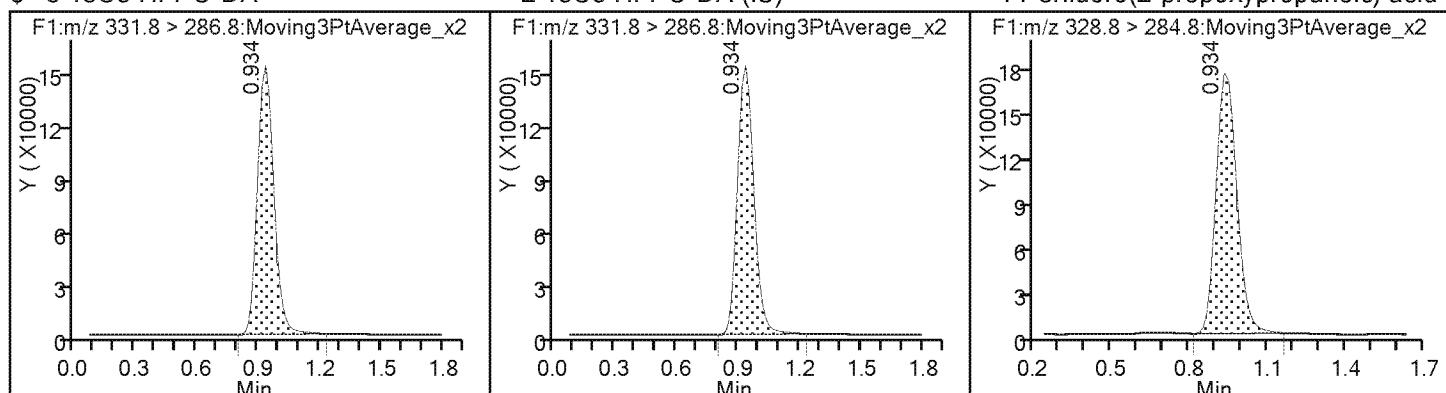

TestAmerica Denver
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30005.d
 Lims ID: LCS 280-403139/2-A
 Client ID:
 Sample Type: LCS
 Inject. Date: 30-Jan-2018 09:19:02 ALS Bottle#: 12 Worklist Smp#: 5
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: LCS280-403139/2-A
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 30-Jan-2018 14:49:56 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM
 Process Host: XAWRK010

First Level Reviewer: meyera Date: 30-Jan-2018 14:47:19

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	10.7	106.82

FORM I
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: TestAmerica Denver Job No.: 320-35355-1
SDG No.:
Client Sample ID: Lab Sample ID: LCSD 280-403139/3-A
Matrix: Water Lab File ID: hfpo718A30006.d
Analysis Method: 8321A Date Collected:
Extraction Method: 3535 Date Extracted: 01/29/2018 21:37
Sample wt/vol: 250 (mL) Date Analyzed: 01/30/2018 09:22
Con. Extract Vol.: 5 (mL) Dilution Factor: 1
Injection Volume: 20 (uL) GC Column: Synergi Hydro ID:
% Moisture: GPC Cleanup: (Y/N) N
Analysis Batch No.: 403246 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	0.207		0.010	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	105		50-200

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30006.d
 Lims ID: LCSD 280-403139/3-A
 Client ID:
 Sample Type: LCSD
 Inject. Date: 30-Jan-2018 09:22:17 ALS Bottle#: 13 Worklist Smp#: 6
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: LCSD280-403139/3-A
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 30-Jan-2018 14:49:56 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK010

First Level Reviewer: meyera Date: 30-Jan-2018 14:47:26

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA
 331.8 > 286.8 0.934 0.961 -0.027 1.000 769072 10.5 1074
 * 2 13C3 HFPO-DA (IS)
 331.8 > 286.8 0.934 0.961 -0.027 1.000 769072 10.0 1074
 1 Perfluoro(2-propoxypropanoic) acid
 328.8 > 284.8 0.934 0.988 -0.054 1.000 822061 10.3 365

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20180130-66815.b\\hfp0718A30006.d

Injection Date: 30-Jan-2018 09:22:17 Instrument ID: LC_LCMS7

Lims ID: LCSD 280-403139/3-A

Client ID:

Operator ID: JBH ALS Bottle#: 13 Worklist Smp#: 6


Injection Vol: 20.0 ul Dil. Factor: 1.0000

Method: HFPO Limit Group: LC - 8321A_HFPO_Du

\$ 3 13C3 HFPO-DA

* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

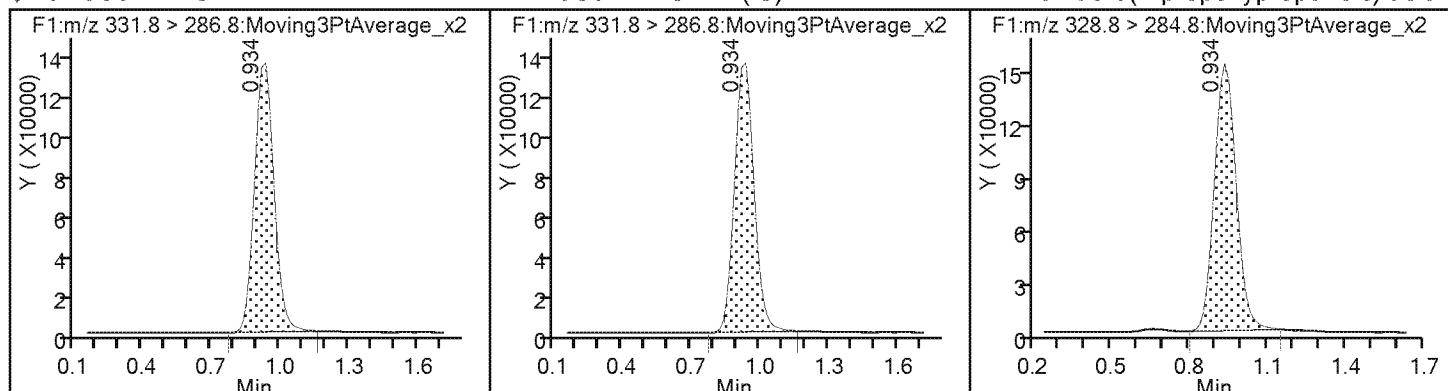

TestAmerica Denver
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30006.d
 Lims ID: LCSD 280-403139/3-A
 Client ID:
 Sample Type: LCSD
 Inject. Date: 30-Jan-2018 09:22:17 ALS Bottle#: 13 Worklist Smp#: 6
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: LCSD280-403139/3-A
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 30-Jan-2018 14:49:56 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK010

First Level Reviewer: meyera Date: 30-Jan-2018 14:47:26

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	10.5	105.14

FORM I
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: TestAmerica Denver Job No.: 320-35355-1

SDG No.: _____

Client Sample ID: _____ Lab Sample ID: LLCS 280-403139/4-A

Matrix: Water Lab File ID: hfpo718A30007.d

Analysis Method: 8321A Date Collected: _____

Extraction Method: 3535 Date Extracted: 01/29/2018 21:37

Sample wt/vol: 250 (mL) Date Analyzed: 01/30/2018 09:25

Con. Extract Vol.: 5 (mL) Dilution Factor: 1

Injection Volume: 20 (uL) GC Column: Synergi Hydro ID: _____

% Moisture: GPC Cleanup: (Y/N) N

Analysis Batch No.: 403246 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	0.0147		0.010	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	128		50-200

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30007.d
 Lims ID: LLCS 280-403139/4-A
 Client ID:
 Sample Type: LLCS
 Inject. Date: 30-Jan-2018 09:25:34 ALS Bottle#: 14 Worklist Smp#: 7
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: LLCS280-403139/4-A
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 30-Jan-2018 14:49:56 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK010

First Level Reviewer: meyera Date: 30-Jan-2018 14:47:40

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA

331.8 > 286.8 0.920 0.961 -0.041 1.000 939514 12.8 1182

* 2 13C3 HFPO-DA (IS)

331.8 > 286.8 0.920 0.961 -0.041 939514 10.0 1182

1 Perfluoro(2-propoxypropanoic) acid M

328.8 > 284.8 0.920 0.988 -0.068 1.000 90632 0.7373 26.8 M

QC Flag Legend

Review Flags

M - Manually Integrated

Report Date: 30-Jan-2018 14:50:00

Chrom Revision: 2.2 24-Jan-2018 15:37:30

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20180130-66815.b\\hfpo718A30007.d

Injection Date: 30-Jan-2018 09:25:34 Instrument ID: LC_LCMS7

Lims ID: LLCS 280-403139/4-A

Client ID:

Operator ID: JBH ALS Bottle#: 14 Worklist Smp#: 7


Injection Vol: 20.0 ul Dil. Factor: 1.0000

Method: HFPO Limit Group: LC - 8321A_HFPO_Du

\$ 3 13C3 HFPO-DA

* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid (M)

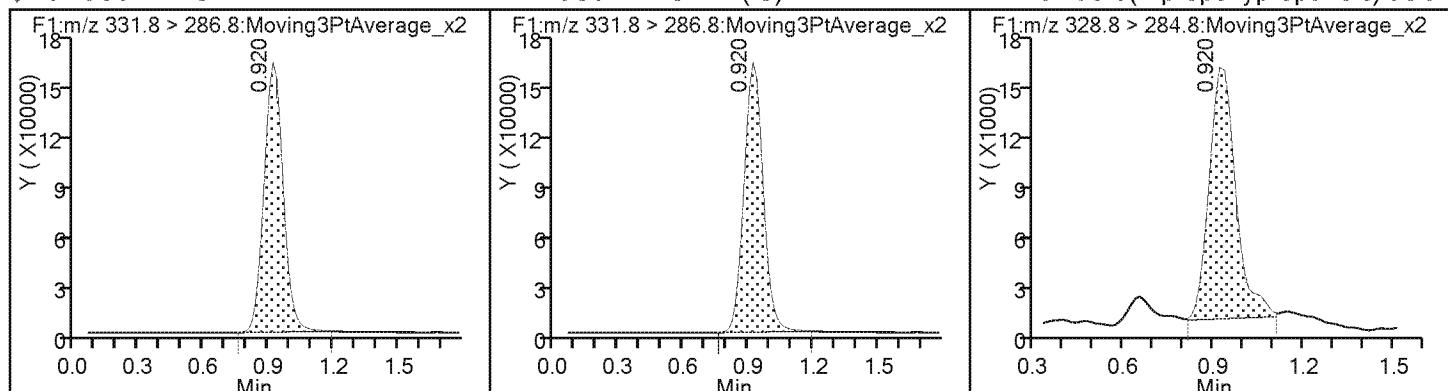

TestAmerica Denver
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30007.d
 Lims ID: LLCS 280-403139/4-A
 Client ID:
 Sample Type: LLCS
 Inject. Date: 30-Jan-2018 09:25:34 ALS Bottle#: 14 Worklist Smp#: 7
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: LLCS280-403139/4-A
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 30-Jan-2018 14:49:56 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM
 Process Host: XAWRK010

First Level Reviewer: meyera Date: 30-Jan-2018 14:47:40

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	12.8	128.45

TestAmerica Denver


Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20180130-66815.b\\hfpo718A30007.d
 Injection Date: 30-Jan-2018 09:25:34 Instrument ID: LC_LCMS7
 Lims ID: LLCS 280-403139/4-A
 Client ID:
 Operator ID: JBH ALS Bottle#: 14 Worklist Smp#: 7
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Method: HFPO Limit Group: LC - 8321A_HFPO_Du
 Column: Detector F1:MRM

1 Perfluoro(2-propoxypropanoic) acid, CAS: 13252-13-6

Signal: 1


RT: 0.92
 Area: 79101
 Amount: 0.616004
 Amount Units: ug/l

Processing Integration Results

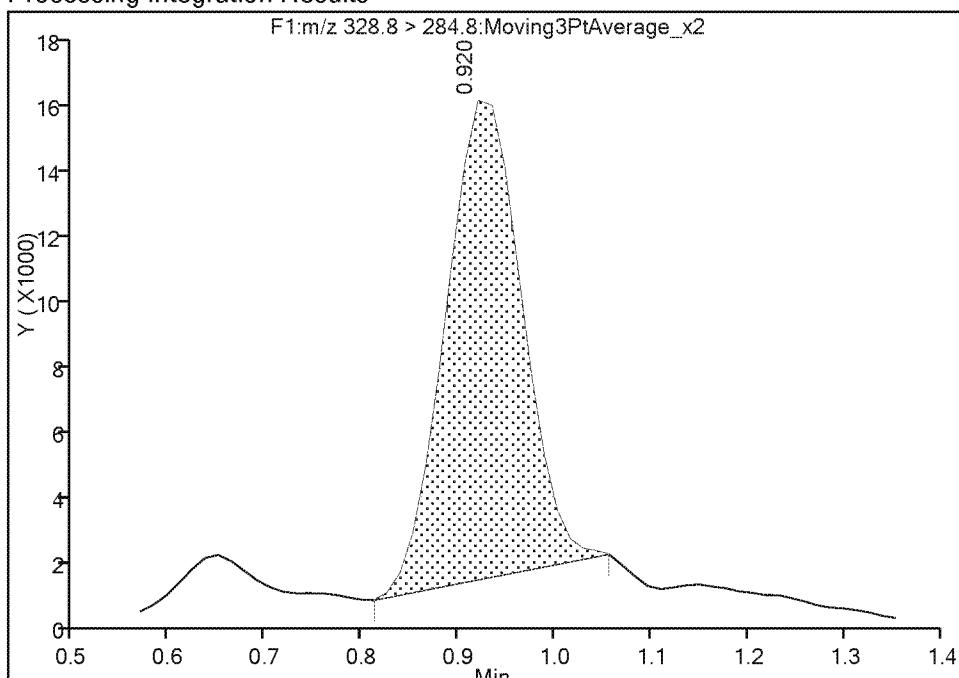

RT: 0.92
 Area: 90632
 Amount: 0.737271
 Amount Units: ug/l

Manual Integration Results

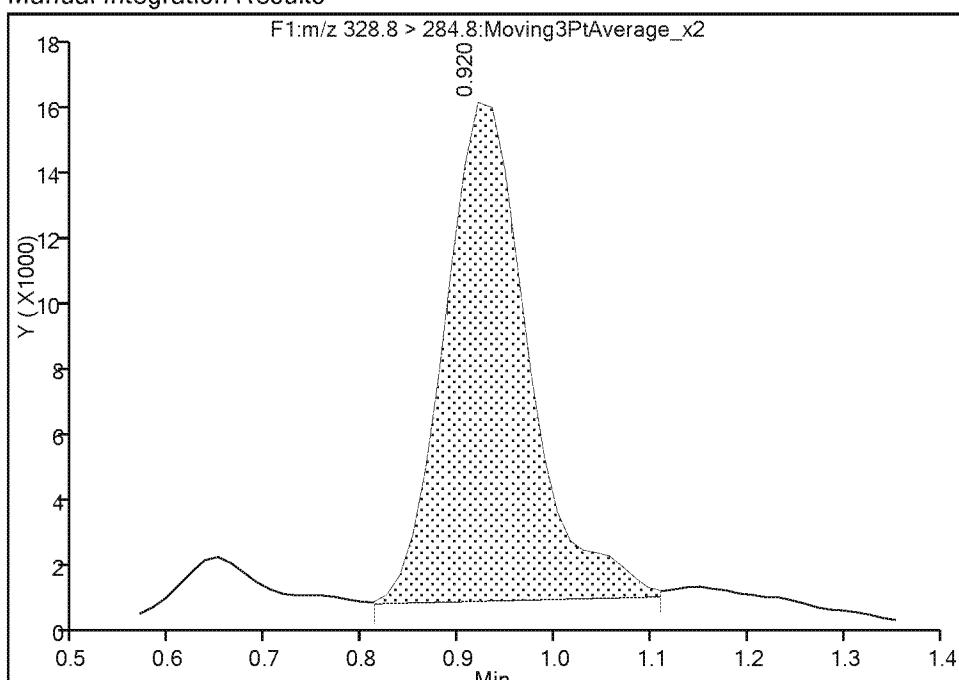

Reviewer: meyera, 30-Jan-2018 14:47:37

Audit Action: Manually Integrated

Audit Reason: Baseline

FORM I
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: TestAmerica Denver Job No.: 320-35355-1

SDG No.: _____

Client Sample ID: _____ Lab Sample ID: DLCK 280-390728/12

Matrix: Water Lab File ID: hfpo717J10035.d

Analysis Method: 8321A Date Collected: _____

Extraction Method: _____ Date Extracted: _____

Sample wt/vol: 1 (mL) Date Analyzed: 10/10/2017 10:04

Con. Extract Vol.: _____ Dilution Factor: 1

Injection Volume: 20 (uL) GC Column: Synergi Hydro ID: _____

% Moisture: _____ GPC Cleanup: (Y/N) N

Analysis Batch No.: 390728 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	<0.50		0.50	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	102		50-200

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10035.d
 Lims ID: DLCK
 Client ID:
 Sample Type: DLCK
 Inject. Date: 10-Oct-2017 10:04:34 ALS Bottle#: 2 Worklist Smp#: 12
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: DLCK
 Misc. Info.: HFPO17J10
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 10-Oct-2017 12:51:51 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK005

First Level Reviewer: meyera Date: 10-Oct-2017 11:51:31

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA
 331.8 > 286.8 0.880 0.880 0.0 1.000 749614 10.2 480
 * 2 13C3 HFPO-DA (IS)
 331.8 > 286.8 0.880 0.880 0.0 1.000 749614 10.0 480
 1 Perfluoro(2-propoxypropanoic) acid
 328.8 > 284.8 0.893 0.885 0.008 1.000 31104 0.1941 16.6

Reagents:

HFPO_CAL-1_00031 Amount Added: 1.00 Units: mL

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20171010-63483.b\\hfp0717J10035.d

Injection Date: 10-Oct-2017 10:04:34 Instrument ID: LC_LCMS7

Lims ID: DLCK

Client ID:

Operator ID: JBH

ALS Bottle#: 2 Worklist Smp#: 12

Injection Vol: 20.0 ul

Dil. Factor: 1.0000


Method: HFPO

Limit Group: LC - 8321A_HFPO_Du

\$ 3 13C3 HFPO-DA

* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

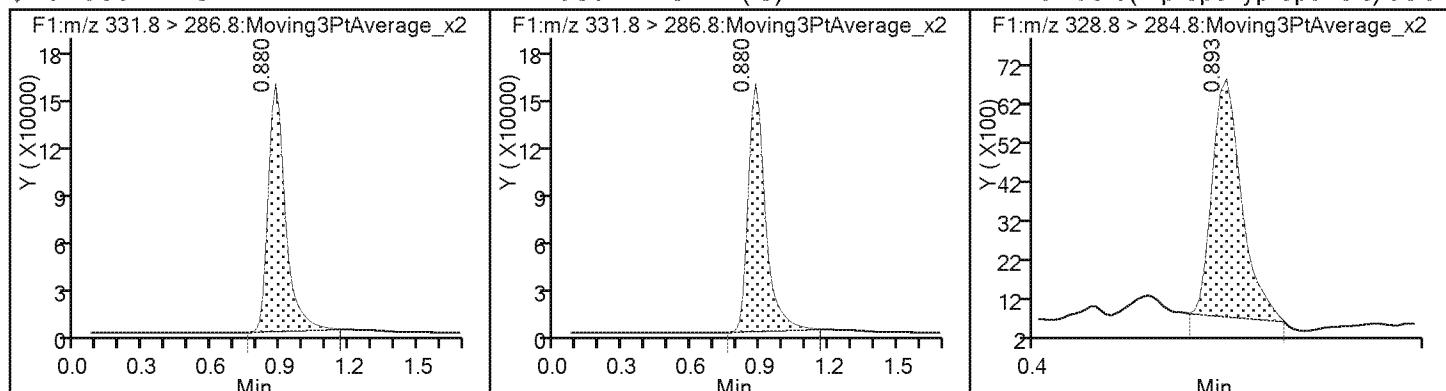

TestAmerica Denver
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10035.d
 Lims ID: DLCK
 Client ID:
 Sample Type: DLCK
 Inject. Date: 10-Oct-2017 10:04:34 ALS Bottle#: 2 Worklist Smp#: 12
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: DLCK
 Misc. Info.: HFPO17J10
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 10-Oct-2017 12:51:51 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM
 Process Host: XAWRK005

First Level Reviewer: meyera Date: 10-Oct-2017 11:51:31

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	10.2	102.48

FORM I
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: TestAmerica Denver Job No.: 320-35355-1

SDG No.: _____

Client Sample ID: FAY-D-7166BUTLE-W1-1-0124 Lab Sample ID: 320-35355-9 MS
18 MS

Matrix: Water Lab File ID: hfpo718A30019.d

Analysis Method: 8321A Date Collected: 01/24/2018 15:35

Extraction Method: 3535 Date Extracted: 01/29/2018 21:37

Sample wt/vol: 246.7 (mL) Date Analyzed: 01/30/2018 10:04

Con. Extract Vol.: 5 (mL) Dilution Factor: 1

Injection Volume: 20 (uL) GC Column: Synergi Hydro ID: _____

% Moisture: GPC Cleanup: (Y/N) N

Analysis Batch No.: 403246 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	0.202		0.010	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	123		50-200

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30019.d
 Lims ID: 320-35355-A-9-A MS
 Client ID: FAY-D-7166BUTLE-W1-1-012418
 Sample Type: MS
 Inject. Date: 30-Jan-2018 10:04:33 ALS Bottle#: 25 Worklist Smp#: 19
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: 320-35355-A-9-AMS
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 05-Feb-2018 13:15:57 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK025

First Level Reviewer: meyera Date: 30-Jan-2018 14:49:07

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA
 331.8 > 286.8 0.907 0.961 -0.054 1.000 898882 12.3 1163
 * 2 13C3 HFPO-DA (IS)
 331.8 > 286.8 0.907 0.961 -0.054 1.000 898882 10.0 1163
 1 Perfluoro(2-propoxypropanoic) acid
 328.8 > 284.8 0.920 0.988 -0.068 1.000 927349 9.98 184

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20180130-66815.b\\hfpo718A30019.d

Injection Date: 30-Jan-2018 10:04:33 Instrument ID: LC_LCMS7

Lims ID: 320-35355-A-9-A MS

Client ID: FAY-D-7166BUTLE-W1-1-012418

Operator ID: JBH ALS Bottle#: 25 Worklist Smp#: 19


Injection Vol: 20.0 ul Dil. Factor: 1.0000

Method: HFPO Limit Group: LC - 8321A_HFPO_Du

\$ 3 13C3 HFPO-DA

* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid

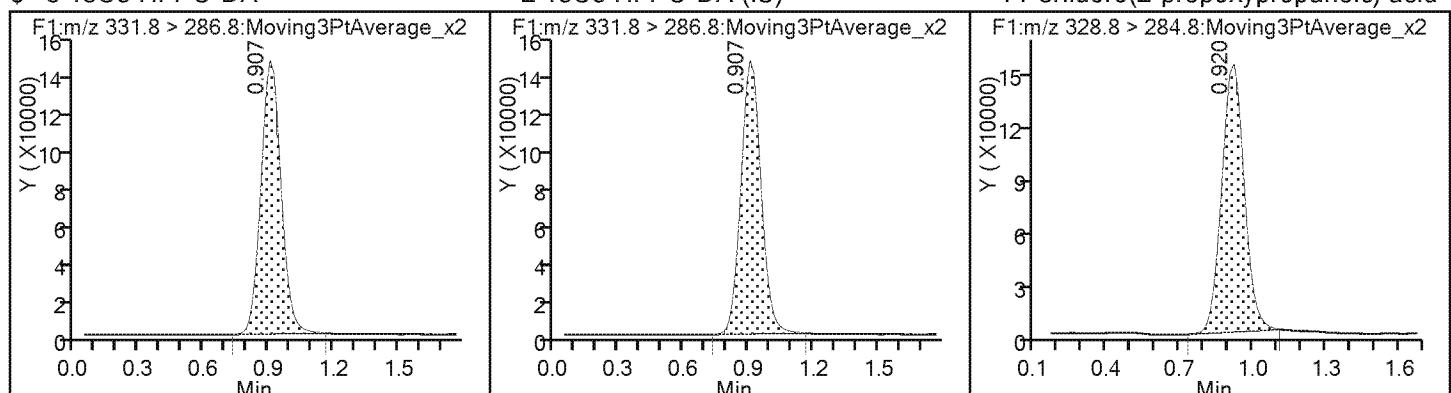

TestAmerica Denver
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30019.d
 Lims ID: 320-35355-A-9-A MS
 Client ID: FAY-D-7166BUTLE-W1-1-012418
 Sample Type: MS
 Inject. Date: 30-Jan-2018 10:04:33 ALS Bottle#: 25 Worklist Smp#: 19
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: 320-35355-A-9-AMS
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 05-Feb-2018 13:15:57 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK025

First Level Reviewer: meyera Date: 30-Jan-2018 14:49:07

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	12.3	122.89

FORM I
LCMS ORGANICS ANALYSIS DATA SHEET

Lab Name: TestAmerica Denver Job No.: 320-35355-1

SDG No.: _____

Client Sample ID: FAY-D-7166BUTLE-W1-1-0124 Lab Sample ID: 320-35355-9 DU
18 DU

Matrix: Water Lab File ID: hfpo718A30018.d

Analysis Method: 8321A Date Collected: 01/24/2018 15:35

Extraction Method: 3535 Date Extracted: 01/29/2018 21:37

Sample wt/vol: 249.1 (mL) Date Analyzed: 01/30/2018 10:01

Con. Extract Vol.: 5 (mL) Dilution Factor: 1

Injection Volume: 20 (uL) GC Column: Synergi Hydro ID: _____

% Moisture: GPC Cleanup: (Y/N) N

Analysis Batch No.: 403246 Units: ug/L

CAS NO.	COMPOUND NAME	RESULT	Q	RL	
13252-13-6	HFPO-DA	0.0171		0.010	

CAS NO.	SURROGATE	%REC	Q	LIMITS
STL02255	13C3 HFPO-DA	100		50-200

TestAmerica Denver
Target Compound Quantitation Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30018.d
 Lims ID: 320-35355-C-9-A DU
 Client ID: FAY-D-7166BUTLE-W1-1-012418
 Sample Type: DU
 Inject. Date: 30-Jan-2018 10:01:18 ALS Bottle#: 24 Worklist Smp#: 18
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: 320-35355-C-9-ADU
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 05-Feb-2018 13:15:57 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK025

First Level Reviewer: meyera Date: 30-Jan-2018 14:49:03

Signal	RT	EXP RT	DLT RT	REL RT	Response	Amount ug/l	S/N	Flags
--------	----	--------	--------	--------	----------	-------------	-----	-------

\$ 3 13C3 HFPO-DA
 331.8 > 286.8 0.907 0.961 -0.054 1.000 733264 10.0 1037
 * 2 13C3 HFPO-DA (IS)
 331.8 > 286.8 0.907 0.961 -0.054 1.000 733264 10.0 1037
 1 Perfluoro(2-propoxypropanoic) acid M
 328.8 > 284.8 0.920 0.988 -0.068 1.000 79179 0.8510 14.2 M

QC Flag Legend

Review Flags

M - Manually Integrated

TestAmerica Denver

Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20180130-66815.b\\hfp0718A30018.d

Injection Date: 30-Jan-2018 10:01:18 Instrument ID: LC_LCMS7

Lims ID: 320-35355-C-9-A DU

Client ID: FAY-D-7166BUTLE-W1-1-012418

Operator ID: JBH ALS Bottle#: 24 Worklist Smp#: 18


Injection Vol: 20.0 ul Dil. Factor: 1.0000

Method: HFPO Limit Group: LC - 8321A_HFPO_Du

\$ 3 13C3 HFPO-DA

* 2 13C3 HFPO-DA (IS)

1 Perfluoro(2-propoxypropanoic) acid (M)

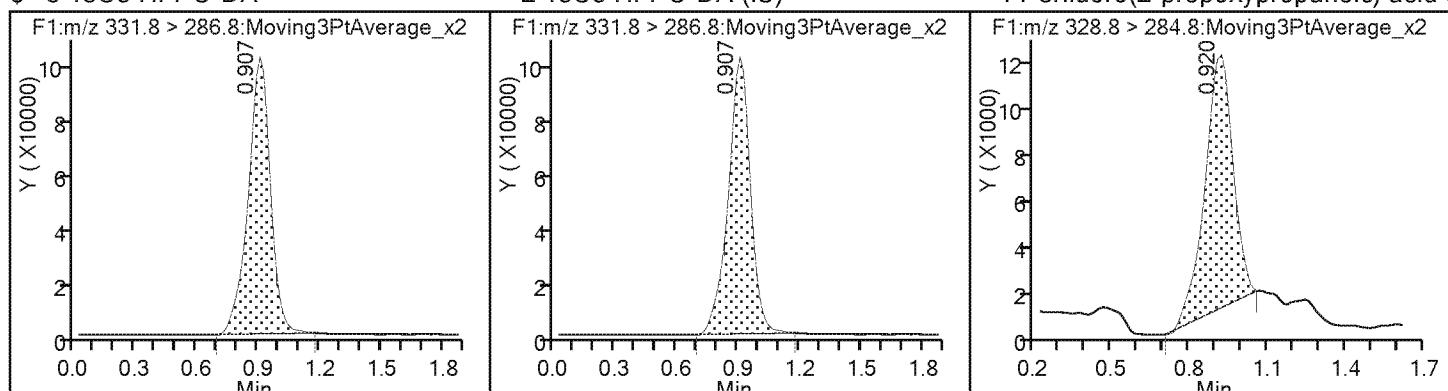

TestAmerica Denver
Recovery Report

Data File: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\hfpo718A30018.d
 Lims ID: 320-35355-C-9-A DU
 Client ID: FAY-D-7166BUTLE-W1-1-012418
 Sample Type: DU
 Inject. Date: 30-Jan-2018 10:01:18 ALS Bottle#: 24 Worklist Smp#: 18
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Sample Info: 320-35355-C-9-ADU
 Misc. Info.: HFPO18A30
 Operator ID: JBH Instrument ID: LC_LCMS7
 Method: \\ChromNA\Denver\ChromData\LC_LCMS7\20180130-66815.b\HFPO.m
 Limit Group: LC - 8321A_HFPO_Du
 Last Update: 05-Feb-2018 13:15:57 Calib Date: 10-Oct-2017 09:58:07
 Integrator: Picker
 Quant Method: Internal/External Standard Quant By: Initial Calibration
 Last ICal File: \\ChromNA\Denver\ChromData\LC_LCMS7\20171010-63483.b\hfpo717J10033.d

Column 1 : Det: F1:MRM

Process Host: XAWRK025

First Level Reviewer: meyera Date: 30-Jan-2018 14:49:03

Compound	Amount Added	Amount Recovered	% Rec.
\$ 3 13C3 HFPO-DA	10.0	10.0	100.25

TestAmerica Denver


Data File: \\ChromNA\\Denver\\ChromData\\LC_LCMS7\\20180130-66815.b\\hfpo718A30018.d
 Injection Date: 30-Jan-2018 10:01:18 Instrument ID: LC_LCMS7
 Lims ID: 320-35355-C-9-A DU
 Client ID: FAY-D-7166BUTLE-W1-1-012418
 Operator ID: JBH ALS Bottle#: 24 Worklist Smp#: 18
 Injection Vol: 20.0 ul Dil. Factor: 1.0000
 Method: HFPO Limit Group: LC - 8321A_HFPO_Du
 Column: Detector F1:MRM

1 Perfluoro(2-propoxypropanoic) acid, CAS: 13252-13-6

Signal: 1


RT: 0.92
 Area: 73336
 Amount: 0.772309
 Amount Units: ug/l

Processing Integration Results

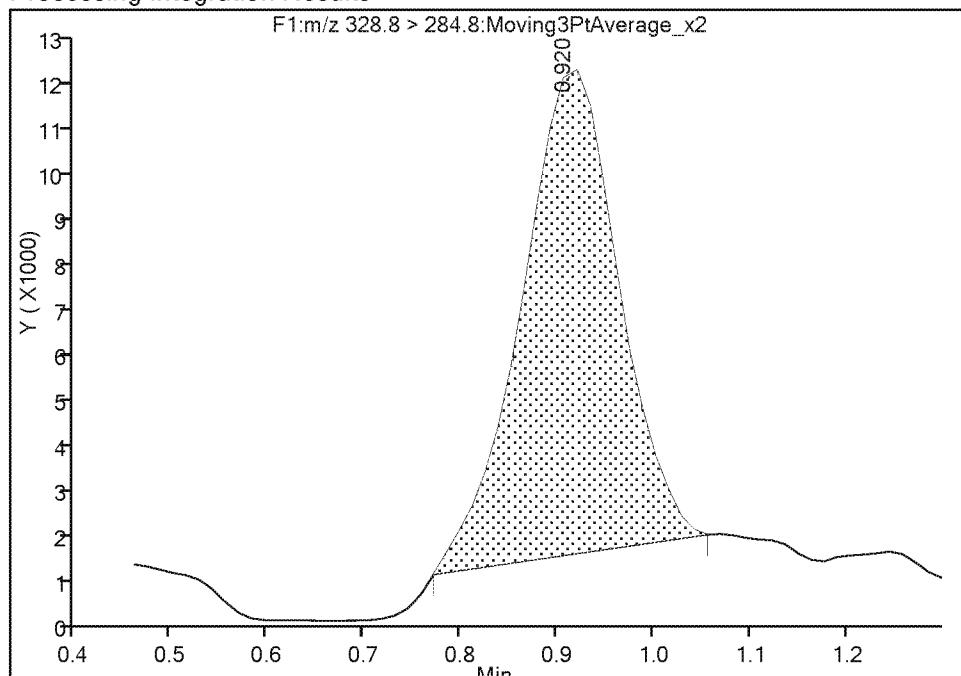

RT: 0.92
 Area: 79179
 Amount: 0.851041
 Amount Units: ug/l

Manual Integration Results

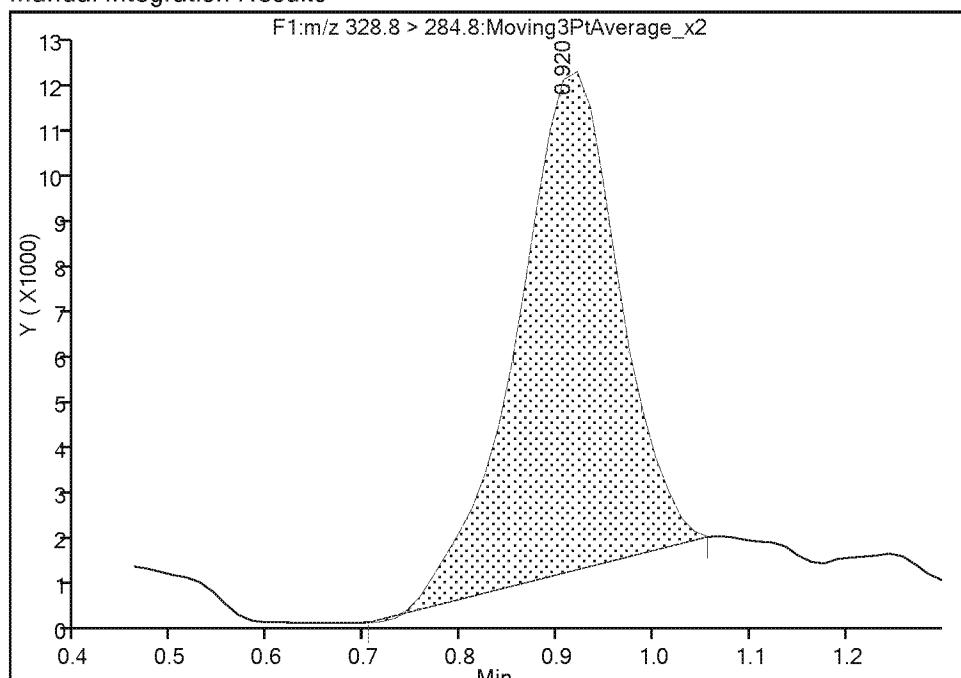

Reviewer: meyera, 05-Feb-2018 13:15:55

Audit Action: Manually Integrated

Audit Reason: Baseline

LCMS ANALYSIS RUN LOG

Lab Name: TestAmerica Denver

Job No.: 320-35355-1

SDG No.:

Instrument ID: LC_LCMS7

Start Date: 09/14/2017 14:40

Analysis Batch Number: 387775

End Date: 09/14/2017 16:05

LAB SAMPLE ID	CLIENT SAMPLE ID	DATE ANALYZED	DILUTION FACTOR	LAB FILE ID	COLUMN ID
STD001 280-387775/3 IC		09/14/2017 14:40	1	hfpo717I14052.d	Synergi Hydro
STD002 280-387775/4 IC		09/14/2017 14:43	1	hfpo717I14053.d	Synergi Hydro
STD003 280-387775/5 IC		09/14/2017 14:46	1	hfpo717I14054.d	Synergi Hydro
STD004 280-387775/6 IC		09/14/2017 14:49	1	hfpo717I14055.d	Synergi Hydro
STD005 280-387775/7 IC		09/14/2017 14:52	1	hfpo717I14056.d	Synergi Hydro
STD006 280-387775/8 IC		09/14/2017 14:55	1	hfpo717I14057.d	Synergi Hydro
STD007 280-387775/9 IC		09/14/2017 14:58	1	hfpo717I14058.d	Synergi Hydro
STD008 280-387775/10 IC		09/14/2017 15:01	1	hfpo717I14059.d	Synergi Hydro
ICB 280-387775/11		09/14/2017 15:04	1		Synergi Hydro
ZZZZZ		09/14/2017 15:07	1		Synergi Hydro
ICV 280-387775/13		09/14/2017 15:10	1	hfpo717I14062.d	Synergi Hydro
ZZZZZ		09/14/2017 15:13	1		Synergi Hydro
ZZZZZ		09/14/2017 15:16	1		Synergi Hydro
ZZZZZ		09/14/2017 15:19	1		Synergi Hydro
ZZZZZ		09/14/2017 15:22	1		Synergi Hydro
ZZZZZ		09/14/2017 15:25	2		Synergi Hydro
ZZZZZ		09/14/2017 15:28	1		Synergi Hydro
ZZZZZ		09/14/2017 15:31	1		Synergi Hydro
ZZZZZ		09/14/2017 15:34	1		Synergi Hydro
ZZZZZ		09/14/2017 15:38	1		Synergi Hydro
ZZZZZ		09/14/2017 15:41	1		Synergi Hydro
CCV 280-387775/24		09/14/2017 15:44	1		Synergi Hydro
ZZZZZ		09/14/2017 15:47	2		Synergi Hydro
ZZZZZ		09/14/2017 15:50	1		Synergi Hydro
ZZZZZ		09/14/2017 15:53	1		Synergi Hydro
ZZZZZ		09/14/2017 15:56	1		Synergi Hydro
ZZZZZ		09/14/2017 15:59	1		Synergi Hydro
ZZZZZ		09/14/2017 16:02	1		Synergi Hydro
CCV 280-387775/31		09/14/2017 16:05	1		Synergi Hydro

LCMS ANALYSIS RUN LOG

Lab Name: TestAmerica Denver

Job No.: 320-35355-1

SDG No.:

Instrument ID: LC_LCMS7

Start Date: 10/10/2017 09:35

Analysis Batch Number: 390728

End Date: 10/10/2017 11:19

LAB SAMPLE ID	CLIENT SAMPLE ID	DATE ANALYZED	DILUTION FACTOR	LAB FILE ID	COLUMN ID
STD001 280-390728/3 IC		10/10/2017 09:35	1	hfpo717J10026.d	Synergi Hydro
STD002 280-390728/4 IC		10/10/2017 09:38	1	hfpo717J10027.d	Synergi Hydro
STD003 280-390728/5 IC		10/10/2017 09:41	1	hfpo717J10028.d	Synergi Hydro
STD004 280-390728/6 IC		10/10/2017 09:45	1	hfpo717J10029.d	Synergi Hydro
STD005 280-390728/7 IC		10/10/2017 09:48	1	hfpo717J10030.d	Synergi Hydro
STD006 280-390728/8 IC		10/10/2017 09:51	1	hfpo717J10031.d	Synergi Hydro
STD007 280-390728/9 IC		10/10/2017 09:54	1	hfpo717J10032.d	Synergi Hydro
STD008 280-390728/10 IC		10/10/2017 09:58	1	hfpo717J10033.d	Synergi Hydro
ICB 280-390728/11		10/10/2017 10:01	1	hfpo717J10034.d	Synergi Hydro
DLCK 280-390728/12		10/10/2017 10:04	1	hfpo717J10035.d	Synergi Hydro
ICV 280-390728/13		10/10/2017 10:07	1	hfpo717J10036.d	Synergi Hydro
ZZZZZ		10/10/2017 10:11	1		Synergi Hydro
ZZZZZ		10/10/2017 10:14	1		Synergi Hydro
ZZZZZ		10/10/2017 10:17	1		Synergi Hydro
ZZZZZ		10/10/2017 10:20	1		Synergi Hydro
ZZZZZ		10/10/2017 10:23	1		Synergi Hydro
ZZZZZ		10/10/2017 10:27	1		Synergi Hydro
ZZZZZ		10/10/2017 10:30	1		Synergi Hydro
ZZZZZ		10/10/2017 10:33	1		Synergi Hydro
ZZZZZ		10/10/2017 10:36	1		Synergi Hydro
ZZZZZ		10/10/2017 10:40	1		Synergi Hydro
CCV 280-390728/24		10/10/2017 10:43	1	hfpo717J10047.d	Synergi Hydro
ZZZZZ		10/10/2017 10:46	1		Synergi Hydro
ZZZZZ		10/10/2017 10:49	1		Synergi Hydro
ZZZZZ		10/10/2017 10:53	1		Synergi Hydro
ZZZZZ		10/10/2017 10:56	1		Synergi Hydro
ZZZZZ		10/10/2017 10:59	1		Synergi Hydro
ZZZZZ		10/10/2017 11:02	1		Synergi Hydro
ZZZZZ		10/10/2017 11:06	1		Synergi Hydro
ZZZZZ		10/10/2017 11:09	1		Synergi Hydro
ZZZZZ		10/10/2017 11:12	1		Synergi Hydro
ZZZZZ		10/10/2017 11:16	1		Synergi Hydro
CCV 280-390728/35		10/10/2017 11:19	1		Synergi Hydro

8321A

LCMS ANALYSIS RUN LOG

Lab Name: TestAmerica Denver

Job No.: 320-35355-1

SDG No.:

Instrument ID: LC_LCMS7

Start Date: 01/30/2018 09:12

Analysis Batch Number: 403246

End Date: 01/30/2018 10:37

LAB SAMPLE ID	CLIENT SAMPLE ID	DATE ANALYZED	DILUTION FACTOR	LAB FILE ID	COLUMN ID
CCV 280-403246/3		01/30/2018 09:12	1	hfpo718A30003.d	Synergi Hydro
MB 280-403139/1-A		01/30/2018 09:15	1	hfpo718A30004.d	Synergi Hydro
LCS 280-403139/2-A		01/30/2018 09:19	1	hfpo718A30005.d	Synergi Hydro
LCSD 280-403139/3-A		01/30/2018 09:22	1	hfpo718A30006.d	Synergi Hydro
LLCS 280-403139/4-A		01/30/2018 09:25	1	hfpo718A30007.d	Synergi Hydro
320-35355-1		01/30/2018 09:28	1	hfpo718A30008.d	Synergi Hydro
320-35355-2		01/30/2018 09:32	1	hfpo718A30009.d	Synergi Hydro
320-35355-3		01/30/2018 09:35	1	hfpo718A30010.d	Synergi Hydro
320-35355-4		01/30/2018 09:38	1	hfpo718A30011.d	Synergi Hydro
CCV 280-403246/12		01/30/2018 09:41	1	hfpo718A30012.d	Synergi Hydro
320-35355-5		01/30/2018 09:45	1	hfpo718A30013.d	Synergi Hydro
320-35355-6		01/30/2018 09:48	1	hfpo718A30014.d	Synergi Hydro
320-35355-7		01/30/2018 09:51	1	hfpo718A30015.d	Synergi Hydro
320-35355-8		01/30/2018 09:54	1	hfpo718A30016.d	Synergi Hydro
320-35355-9		01/30/2018 09:58	1	hfpo718A30017.d	Synergi Hydro
320-35355-9 DU		01/30/2018 10:01	1	hfpo718A30018.d	Synergi Hydro
320-35355-9 MS		01/30/2018 10:04	1	hfpo718A30019.d	Synergi Hydro
320-35355-10		01/30/2018 10:07	1	hfpo718A30020.d	Synergi Hydro
CCV 280-403246/21		01/30/2018 10:11	1	hfpo718A30021.d	Synergi Hydro
ZZZZZ		01/30/2018 10:14	1		Synergi Hydro
ZZZZZ		01/30/2018 10:17	1		Synergi Hydro
ZZZZZ		01/30/2018 10:20	1		Synergi Hydro
ZZZZZ		01/30/2018 10:24	1		Synergi Hydro
ZZZZZ		01/30/2018 10:27	1		Synergi Hydro
ZZZZZ		01/30/2018 10:30	1		Synergi Hydro
ZZZZZ		01/30/2018 10:33	1		Synergi Hydro
CCV 280-403246/29		01/30/2018 10:37	1		Synergi Hydro

LCMS BATCH WORKSHEET

Lab Name: TestAmerica Denver

Job No.: 320-35355-1

SDG No.:

Batch Number: 403139

Batch Start Date: 01/29/18 21:37

Batch Analyst: Cokley, Cheyana D

Batch Method: 3535

Batch End Date: 01/30/18 01:54

Lab Sample ID	Client Sample ID	Method Chain	Basis	GrossWeight	TareWeight	InitialAmount	FinalAmount	HFPO I.S. 00007	HFPO Spike 00004
MB 280-403139/1		3535, 8321A				250 mL	5 mL	0.1 mL	
LCS 280-403139/2		3535, 8321A				250 mL	5 mL	0.1 mL	0.1 mL
LCSD 280-403139/3		3535, 8321A				250 mL	5 mL	0.1 mL	0.1 mL
LLCS 280-403139/4		3535, 8321A				250 mL	5 mL	0.1 mL	0.01 mL
320-35355-B-1	FAY-D-7154BUTLE-W1-1-012418	3535, 8321A	T	279.0 g	27.4 g	251.6 mL	5 mL	0.1 mL	
320-35355-B-2	FAY-D-384JERML-W1-1-012418	3535, 8321A	T	280.0 g	27.0 g	253 mL	5 mL	0.1 mL	
320-35355-B-3	FAY-D-2528INDIA-W1-2-012418	3535, 8321A	T	278.5 g	27.2 g	251.3 mL	5 mL	0.1 mL	
320-35355-B-4	FAY-D-661PINEA-W1-1-012418	3535, 8321A	T	264.1 g	27.7 g	236.4 mL	5 mL	0.1 mL	
320-35355-B-5	FAY-D-701PINEA-W1-1-012418	3535, 8321A	T	268.2 g	27.7 g	240.5 mL	5 mL	0.1 mL	
320-35355-B-6	FAY-D-701PINEA-W1-2-012418	3535, 8321A	T	279.3 g	27.0 g	252.3 mL	5 mL	0.1 mL	
320-35355-B-7	FAY-D-663PINEA-W1-1-012418	3535, 8321A	T	270.0 g	27.9 g	242.1 mL	5 mL	0.1 mL	
320-35355-A-8	FAY-D-7166BUTLE-W2-1-012418	3535, 8321A	T	283.6 g	27.5 g	256.1 mL	5 mL	0.1 mL	
320-35355-B-9	FAY-D-7166BUTLE-W1-1-012418	3535, 8321A	T	280.8 g	28.5 g	252.3 mL	5 mL	0.1 mL	
320-35355-C-9	FAY-D-7166BUTLE-W1-1-012418	3535, 8321A	T	278.0 g	28.9 g	249.1 mL	5 mL	0.1 mL	
320-35355-A-9	FAY-D-7166BUTLE-W1-1-012418	3535, 8321A	T	274.4 g	27.7 g	246.7 mL	5 mL	0.1 mL	0.1 mL
320-35355-D-10	FAY-D-7166BUTLE-W1-1-012418D	3535, 8321A	T	283.2 g	26.9 g	256.3 mL	5 mL	0.1 mL	

The pound sign (#) in the amount added field denotes that the reagent was used undiluted. All calculations are performed using the stated concentration for this reagent.

LCMS BATCH WORKSHEET

Lab Name: TestAmerica Denver

Job No.: 320-35355-1

SDG No.:

Batch Number: 403139

Batch Start Date: 01/29/18 21:37

Batch Analyst: Cokley, Cheyana D

Batch Method: 3535

Batch End Date: 01/30/18 01:54

Batch Notes	
Acid ID	2% Formic Aci_00140
Acid Name	2% Formic Acid
Balance ID	24350888
Batch Comment	Reviewer:cdc
First End time	1.29@2250 lines 21-23 on 1.30.18@0045
H2O ID	HPLC_Water_00847
Pipette ID	F, SPE-1+ syringe
Reagent ID	10% NH4OH
Reagent Lot Number	10% NH4OH_00117
Solvent Lot #	HPLC_MeOH_00016
Solvent Name	Methanol
SOP Number	DV-OP-0019
SPE Cartridge Type	STRATA-X-AW (8B S038 FCH)
Solid Phase Extraction Disk ID	S308-0076
First Start time	1.29@2210 lines 21-23 @2358

Basis	Basis Description
T	Total/NA

The pound sign (#) in the amount added field denotes that the reagent was used undiluted. All calculations are performed using the stated concentration for this reagent.



Reagent ID: **HFPO_CAL-6_00079**

Description:	level6	Expiration Date:	02/09/2018
No. of Bottles:	1	Laboratory:	TestAmerica Denver
Storage Location:	LCMS	Prepared By:	Meyer, Andrew GC
Reagent Volume:	1.000 mL	Solvent:	80:20 Methanol : H2O
Creation Date:	01/26/2018	Solvent Lot:	00018
Open Date:			
Container(s):	4938570		
Comment:	level-6		

Reagent Analyte Information

Analyte	Source ID	Source Exp. Date	Source Conc.	Source Conc. Units	Final Conc.	Final Conc. Units
13C3 HFPO-DA	HFPO.I.S._00004	03/28/2018	0.60000	ug/mL	10.00000	ug/L
13C3 HFPO-DA (S)	HFPO.I.S._00004	03/28/2018	0.60000	ug/mL	10.00000	ug/L
Perfluoro(2-propoxypropanoic) acid	HFPO_Spike_00004	10/30/2018	0.60000	ug/mL	6.00000	ug/L

Source Reagents

Reagent	Description	Type	Expiration	Vendor	Vendor Lot #	Vendor Cat Lot #	Volume Used	Volume Units
HFPO.I.S._00004	Internal Standard for HFPO 0.6ug/ml		03/28/18				20.00000	uL
HFPO_Spike_00004	HFPO LC/MS Calibration Spike 0.6ug/ml		10/30/18				10.00000	uL



Reagent ID: **HFPO_CAL-6_00079**

Description:	level6	Expiration Date:	02/09/2018
No. of Bottles:	1	Laboratory:	TestAmerica Denver
Storage Location:	LCMS	Prepared By:	Meyer, Andrew GC
Reagent Volume:	1.000 mL	Solvent:	60:20 Methanol : H2O
Creation Date:	01/28/2018	Solvent Lot:	00018
Open Date:			
Container(s):	4938671		
Comment:	level6		

Reagent Analyte Information

Analyte	Source ID	Source Exp. Date	Source Conc.	Source Conc. Units	Final Conc.	Final Conc. Units
13C3 HFPO-DA	HFPO I.S._00004	09/28/2018	0.50000	ug/mL	10.00000	ug/L
13C3 HFPO-DA (IS)	HFPO I.S._00004	09/28/2018	0.50000	ug/mL	10.00000	ug/L
Perfluoro(2-propoxypropanoic) acid	HFPO Spike_00004	10/30/2018	0.50000	ug/mL	10.00000	ug/L

Source Reagents

Reagent	Description	Type	Expiration	Vendor	Vendor Lot #	Vendor Cat Lot #	Volume Used	Volume Units
HFPO I.S._00004	Internal Standard for HFPO 0.5ug/ml		09/28/18				20.00000	uL
HFPO Spike_00004	HFPO LC/Calibration Spike 0.5ug/ml		10/30/18				20.00000	uL

Shipping and Receiving Documents

Chain of Custody Record

ED 002011 00000054-00226

Login Sample Receipt Checklist

Client: Chemours Company FC, LLC The

Job Number: 320-35355-1

Login Number: 35355

List Source: TestAmerica Sacramento

List Number: 1

Creator: Turpen, Troy

Question	Answer	Comment
Radioactivity wasn't checked or is </= background as measured by a survey meter.	N/A	
The cooler's custody seal, if present, is intact.	True	
Sample custody seals, if present, are intact.	N/A	
The cooler or samples do not appear to have been compromised or tampered with.	True	
Samples were received on ice.	True	
Cooler Temperature is acceptable.	True	
Cooler Temperature is recorded.	True	
COC is present.	True	
COC is filled out in ink and legible.	True	
COC is filled out with all pertinent information.	True	
Is the Field Sampler's name present on COC?	True	
There are no discrepancies between the containers received and the COC.	True	
Samples are received within Holding Time (excluding tests with immediate HTs)	True	
Sample containers have legible labels.	True	
Containers are not broken or leaking.	True	
Sample collection date/times are provided.	True	
Appropriate sample containers are used.	True	
Sample bottles are completely filled.	True	
Sample Preservation Verified.	N/A	
There is sufficient vol. for all requested analyses, incl. any requested MS/MSDs	True	
Containers requiring zero headspace have no headspace or bubble is <6mm (1/4").	N/A	
Multiphasic samples are not present.	True	
Samples do not require splitting or compositing.	True	
Residual Chlorine Checked.	N/A	

Login Sample Receipt Checklist

Client: Chemours Company FC, LLC The

Job Number: 320-35355-1

Login Number: 35355

List Source: TestAmerica Denver

List Number: 2

List Creation: 01/27/18 02:42 PM

Creator: Pottruff, Reed W

Question	Answer	Comment
Radioactivity wasn't checked or is </= background as measured by a survey meter.	N/A	
The cooler's custody seal, if present, is intact.	True	
Sample custody seals, if present, are intact.	True	
The cooler or samples do not appear to have been compromised or tampered with.	True	
Samples were received on ice.	True	
Cooler Temperature is acceptable.	True	
Cooler Temperature is recorded.	True	
COC is present.	True	
COC is filled out in ink and legible.	True	
COC is filled out with all pertinent information.	True	
Is the Field Sampler's name present on COC?	N/A	Received project as a subcontract.
There are no discrepancies between the containers received and the COC.	True	
Samples are received within Holding Time (excluding tests with immediate HTs)	True	
Sample containers have legible labels.	True	
Containers are not broken or leaking.	True	
Sample collection date/times are provided.	True	
Appropriate sample containers are used.	True	
Sample bottles are completely filled.	True	
Sample Preservation Verified.	N/A	
There is sufficient vol. for all requested analyses, incl. any requested MS/MSDs	True	
Containers requiring zero headspace have no headspace or bubble is <6mm (1/4").	N/A	
Multiphasic samples are not present.	True	
Samples do not require splitting or compositing.	True	
Residual Chlorine Checked.	N/A	

Login Sample Receipt Checklist

Client: Chemours Company FC, LLC The

Job Number: 320-35355-1

Login Number: 35355

List Source: TestAmerica Denver

List Number: 3

List Creation: 01/29/18 12:44 PM

Creator: Bieniulis, Dylan T

Question	Answer	Comment
Radioactivity wasn't checked or is </= background as measured by a survey meter.	N/A	
The cooler's custody seal, if present, is intact.	True	
Sample custody seals, if present, are intact.	True	
The cooler or samples do not appear to have been compromised or tampered with.	True	
Samples were received on ice.	True	
Cooler Temperature is acceptable.	True	
Cooler Temperature is recorded.	True	
COC is present.	True	
COC is filled out in ink and legible.	True	
COC is filled out with all pertinent information.	True	
Is the Field Sampler's name present on COC?	True	
There are no discrepancies between the containers received and the COC.	True	
Samples are received within Holding Time (excluding tests with immediate HTs)	True	
Sample containers have legible labels.	True	
Containers are not broken or leaking.	True	
Sample collection date/times are provided.	True	
Appropriate sample containers are used.	True	
Sample bottles are completely filled.	True	
Sample Preservation Verified.	N/A	
There is sufficient vol. for all requested analyses, incl. any requested MS/MSDs	True	
Containers requiring zero headspace have no headspace or bubble is <6mm (1/4").	N/A	
Multiphasic samples are not present.	True	
Samples do not require splitting or compositing.	True	
Residual Chlorine Checked.	N/A	