

From: Strynar, Mark [Strynar.Mark@epa.gov]
Sent: 7/22/2020 8:02:23 PM
To: Linda Birnbaum [Ex. 6 Personal Privacy (PP)]
Subject: Re: Question re to Chemours

Can I get this paper Mottaleb et al., 2019?

Mark

Dr. Mark Strynar
US EPA/ORD/CEMM
strynar.mark@epa.gov
(office) 919-541-3706
(mobile) [Ex. 6 Personal Privacy (PP)]

From: Linda Birnbaum [Ex. 6 Personal Privacy (PP)]
Sent: Wednesday, July 22, 2020 4:01 PM
To: Strynar, Mark <Strynar.Mark@epa.gov>
Subject: RE: Question re to Chemours

Here's the description, so far, in the paper we're working on re to the method....Please don't share yet..

PFAS extraction and quantitation


PFAS were extracted and analyzed as previously described (Mottaleb et al., 2019). Briefly, 50 µL of serum was extracted in 96-well Impact Protein Precipitation Plates (Phenomenex, CA, USA) in methanol with known concentrations of 13C isotope labeled PFAS compounds (PFOSL, PFOAL, PFNAL and FTSL) as surrogate controls. Dried extract was reconstituted in ammonium acetate that contained 3.0 ng/mL of 13C4–PFHpA as internal standard. All extraction plates contained reagent blank (R-BLK), IS-blank and quality control samples. PFAS were quantitated by UHPLC-MS/MS with a Shimadzu ultra-performance liquid chromatography (UPLC; model: Nexra X2 LC 30 AD) coupled to a Sciex QTRAP mass spectrometer (MS; model: Sciex 6500 plus) collecting in negative multiple reaction monitoring (MRM) mode. Parent and product ions for all PFAS, surrogates, and I.S. as well as instrument settings, and analytical method performance parameters have been reported previously (Mottaleb et al., 2019). PFAS were quantitated using IS-based calibration curves with limits of quantitation between 0.04 ng/mL and 0.20 ng/mL depending on the chemical (Mottaleb et al., 2019).

Love your thoughts...

From: Strynar, Mark <Strynar.Mark@epa.gov>
Sent: Wednesday, July 22, 2020 3:55 PM
To: Linda Birnbaum [Ex. 6 Personal Privacy (PP)]
Subject: Re: Question re to Chemours

I can measure what is in it. FYI I found
this <http://adem.alabama.gov/programs/water/drinkingwater/files/Q1PFASDetections.pdf>

Coosa Valley Water Supply District is close to Anniston, AL. Their numbers I would say look suspect with the
PFBS being as high as they are. I will keep digging.
Mark

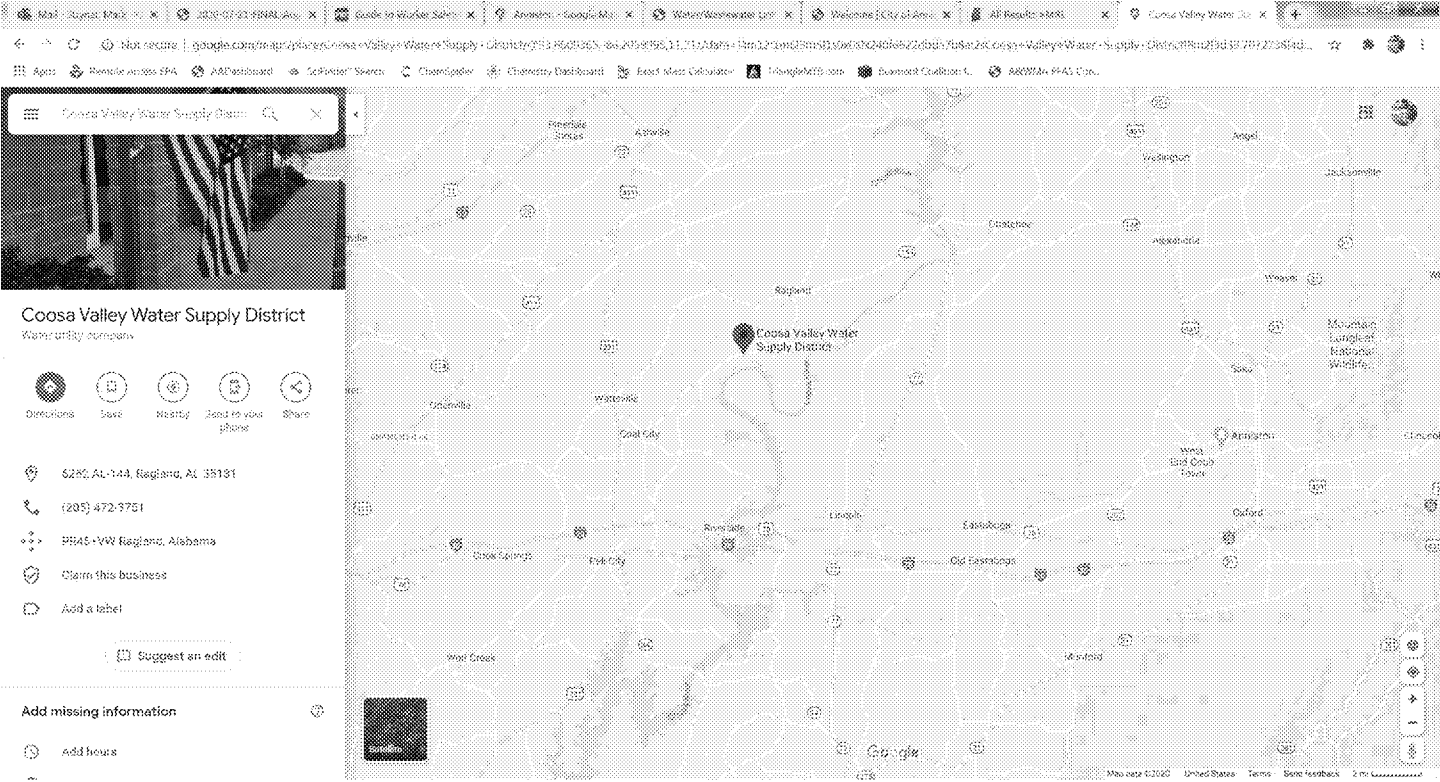

Is this in serum? We have not yet found it in any human serum samples thus far. Could I check out the data and supporting info you are getting with these detection's.

As far as where it may be coming from do you know if the people get drinking water from the Coosa river. I see it flows right by the area. The Coosa starts in Rome GA which has been shown in the past to be contaminated with PFAS due to work going on upstream in Dalton, GA the carpet capital of the world. Can you get a drinking water sample to check for PFAS from the area? It is way upstream but just a thing to rule out.

I don't know if any waste trucks come by this area that may haul HPFO-DA.

Mark

Dr. Mark Strynar
US EPA/ORD/CEMM
strynar.mark@epa.gov
(office) 919-541-3706
(mobile)

Ex. 6 Personal Privacy (PP)

From: Linda Birnbaum

Ex. 6 Personal Privacy (PP)

 >
Sent: Wednesday, July 22, 2020 2:56 PM
To: Strynar, Mark <Strynar.Mark@epa.gov>
Subject: Question re to Chemours

Not sure if I've told you that we have found GenX in ~9% of our cohort in Anniston, AL. Sue Fenton thought you might know if the Chemours truck stps there. We have no idea where these people would get it from – and it was ~2% of population in 2005-7. This is our study looking at PCBs, dioxins, etc. in Anniston – place where PCBs were made from 1929—71....

Linda S. Birnbaum, Ph.D., D.A.B.T., A.T.S.
Scientist Emeritus and Former Director
National Institute of Environmental Health Sciences and National Toxicology Program

Ex. 6 Personal Privacy (PP)

919-280-2884