

CERTIFICATION

SDG No: JC19023 Laboratory: Accutest, New Jersey
Accutest, Florida
Site: BMS, Building 5 Area, PR Matrix: Soil
Humacao, PR

SUMMARY: Soil samples (Table 1) were collected on the BMSMC facility – Building 5 Area. The BMSMC facility is located in Humacao, PR. Samples were taken April 22-25, 2016 and were analyzed in Accutest Laboratory of Dayton, New Jersey for the ABN TCL Special List and for TCL pesticides list that reported the data under SDG No.: JC19023. Accutest Laboratory of Orlando, Florida analyzed for low molecular weight alcohols (LMWA) that also reported the data under SDG No.: JC19023. Results were validated using the latest validation guidelines (July, 2015) of the EPA Hazardous Waste Support Section. The analyses performed are shown in Table 1. Individual data review worksheets are enclosed for each target analyte group. The data sample organic data samples summary form shows for analytes results that were qualified.

In summary the results are valid and can be used for decision taking purposes.

Table 1. Samples analyzed and analysis performed

SAMPLE ID	SAMPLE DESCRIPTION	MATRIX	ANALYSIS PERFORMED
JC19023-1	RA17 (17-18)	Soil	ABN TCL special list; pesticides TCL list; LMWA
JC19023-2	S-40S (10.5 – 11.5)	Soil	ABN TCL special list; pesticides TCL list; LMWA
JC19023-3	BPEB-6	AQ – Equipment Blank	ABN TCL special list; pesticides TCL list; LMWA
JC19023-4	S-39 (11 – 12)	Soil	ABN TCL special list; pesticides TCL list; LMWA

Reviewer Name: Rafael Infante
Chemist License 1888

Signature:

Rafael Infante

Date:

May 17, 2016



SGS Accutest

Report of Analysis

Page 1 of 3

Client Sample ID:	RA-17 (17-18)	Date Sampled:	04/22/16
Lab Sample ID:	JC19023-1	Date Received:	04/26/16
Matrix:	SO - Soil	Percent Solids:	72.6
Method:	SW846 8270D SW846 3546		
Project:	BMSMC, Building 5 Area, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	M124228.D	1	05/06/16	AD	04/29/16	OP93478	EM5253
Run #2							

Run #	Initial Weight	Final Volume
Run #1	31.4 g	1.0 ml
Run #2		

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units	Q
95-57-8	2-Chlorophenol	ND	88	30	ug/kg	
59-50-7	4-Chloro-3-methyl phenol	ND	220	40	ug/kg	
120-83-2	2,4-Dichlorophenol	ND	220	43	ug/kg	
105-67-9	2,4-Dimethylphenol	ND	220	110	ug/kg	
51-28-5	2,4-Dinitrophenol	ND	220	170	ug/kg	
534-52-1	4,6-Dinitro-o-cresol	ND	220	64	ug/kg	
95-48-7	2-Methylphenol	ND	88	50	ug/kg	
	3&4-Methylphenol	ND	88	50	ug/kg	
88-75-5	2-Nitrophenol	ND	220	38	ug/kg	
100-02-7	4-Nitrophenol	ND	440	120	ug/kg	
87-86-5	Pentachlorophenol	ND	220	37	ug/kg	
108-95-2	Phenol	ND	88	34	ug/kg	
58-90-2	2,3,4,6-Tetrachlorophenol	ND	220	57	ug/kg	
95-95-4	2,4,5-Trichlorophenol	ND	220	63	ug/kg	
88-06-2	2,4,6-Trichlorophenol	ND	220	48	ug/kg	
83-32-9	Acenaphthene	ND	44	8.3	ug/kg	
208-96-8	Acenaphthylene	ND	44	6.1	ug/kg	
98-86-2	Acetophenone	ND	220	22	ug/kg	
120-12-7	Anthracene	ND	44	19	ug/kg	
1912-24-9	Atrazine	ND	88	14	ug/kg	
56-55-3	Benzo(a)anthracene	ND	44	6.8	ug/kg	
50-32-8	Benzo(a)pyrene	ND	44	11	ug/kg	
205-99-2	Benzo(b)fluoranthene	ND	44	9.2	ug/kg	
191-24-2	Benzo(g,h,i)perylene	ND	44	12	ug/kg	
207-08-9	Benzo(k)fluoranthene	ND	44	12	ug/kg	
101-55-3	4-Bromophenyl phenyl ether	ND	88	20	ug/kg	
85-68-7	Butyl benzyl phthalate	ND	88	25	ug/kg	
92-52-4	1,1'-Biphenyl	ND	88	12	ug/kg	
100-52-7	Benzaldehyde	103	220	11	ug/kg	
91-58-7	2-Chloronaphthalene	ND	88	8.3	ug/kg	
106-47-8	4-Chloroaniline	ND	220	13	ug/kg	
86-74-8	Carbazole	ND	88	8.6	ug/kg	

ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound



Report of Analysis

Client Sample ID: RA-17 (17-18)
 Lab Sample ID: JC19023-1
 Matrix: SO - Soil
 Method: SW846 8270D SW846 3546
 Project: BSMC, Building 5 Area, PR

Date Sampled: 04/22/16
 Date Received: 04/26/16
 Percent Solids: 72.6

AEN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units	Q
105-60-2	Caprolactam	ND	88	32	ug/kg	
218-01-9	Chrysene	ND	44	11	ug/kg	
111-91-1	bis(2-Chloroethoxy)methane	ND	88	9.3	ug/kg	
111-44-4	bis(2-Chloroethyl)ether	ND	88	25	ug/kg	
108-60-1	bis(2-Chloroisopropyl)ether	ND	88	17	ug/kg	
7005-72-3	4-Chlorophenyl phenyl ether	ND	88	12	ug/kg	
121-14-2	2,4-Dinitrotoluene	ND	44	16	ug/kg	
606-20-2	2,6-Dinitrotoluene	ND	44	13	ug/kg	
91-94-1	3,3'-Dichlorobenzidine	ND	88	59	ug/kg	
53-70-3	Dibenzo(a,h)anthracene	ND	44	17	ug/kg	
132-64-9	Dibenzofuran	ND	88	7.9	ug/kg	
84-74-2	Di-n-butyl phthalate	ND	88	27	ug/kg	
117-84-0	Di-n-octyl phthalate	ND	88	28	ug/kg	
84-66-2	Diethyl phthalate	ND	88	9.9	ug/kg	
131-11-3	Dimethyl phthalate	ND	88	9.0	ug/kg	
117-81-7	bis(2-Ethylhexyl)phthalate	ND	88	9.6	ug/kg	
206-44-0	Fluoranthene	24.3	44	20	ug/kg	J
86-73-7	Fluorene	ND	44	17	ug/kg	
118-74-1	Hexachlorobenzene	ND	88	11	ug/kg	
87-68-3	Hexachlorobutadiene	ND	44	17	ug/kg	
77-47-4	Hexachlorocyclopentadiene	ND	440	26	ug/kg	
67-72-1	Hexachloroethane	ND	220	33	ug/kg	
193-39-5	Indeno(1,2,3-cd)pyrene	ND	44	15	ug/kg	
78-59-1	Isophorone	ND	88	9.7	ug/kg	
90-12-0	1-Methylnaphthalene	ND	88	7.8	ug/kg	
91-57-6	2-Methylnaphthalene	ND	88	36	ug/kg	
88-74-4	2-Nitroaniline	ND	220	32	ug/kg	
99-09-2	3-Nitroaniline	ND	220	15	ug/kg	
100-01-6	4-Nitroaniline	ND	220	17	ug/kg	
98-95-3	Nitrobenzene	ND	88	19	ug/kg	
621-64-7	N-Nitroso-di-n-propylamine	ND	88	20	ug/kg	
86-30-6	N-Nitrosodiphenylamine	ND	220	25	ug/kg	
85-01-8	Phenanthrene	33.9	44	10	ug/kg	J
129-00-0	Pyrene	19.1	44	7.7	ug/kg	J
95-94-3	1,2,4,5-Tetrachlorobenzene	ND	220	9.8	ug/kg	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
367-12-4	2-Fluorophenol	68%		30-106%
4165-62-2	Phenol-d5	71%		30-106%

ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound



Report of Analysis

Client Sample ID:	RA-17 (17-18)	Date Sampled:	04/22/16
Lab Sample ID:	JC19023-1	Date Received:	04/26/16
Matrix:	SO - Soil	Percent Solids:	72.6
Method:	SW846 8270D SW846 3546		
Project:	BMSMC, Building 5 Area, PR		

ABN TCL Special List

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
118-79-6	2,4,6-Tribromophenol	75%		24-140%
4165-60-0	Nitrobenzene-d5	49%		26-122%
321-60-8	2-Fluorobiphenyl	64%		36-112%
1718-51-0	Terphenyl-d14	66%		36-132%



ND = Not detected MDL = Method Detection Limit
RL = Reporting Limit
E = Indicates value exceeds calibration range

J = Indicates an estimated value
B = Indicates analyte found in associated method blank
N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID:	RA-17 (17-18)	Date Sampled:	04/22/16
Lab Sample ID:	JC19023-1	Date Received:	04/26/16
Matrix:	SO - Soil	Percent Solids:	72.6
Method:	SW846 8270D BY SIM SW846 3546		
Project:	BMSMC, Building 5 Area, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	4M65189.D	1	05/09/16	LK	04/29/16	OP93478A	E4M2910
Run #2							

Run #	Initial Weight	Final Volume
Run #1	31.4 g	1.0 ml
Run #2		

CAS No.	Compound	Result	RL	MDL	Units	Q
123-91-1	1,4-Dioxane ^a	7.80	4.4	0.88	ug/kg	
91-20-3	Naphthalene	30.0	4.4	0.54	ug/kg	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
4165-60-0	Nitrobenzene-d5	51%		15-138%
321-60-8	2-Fluorobiphenyl	60%		12-148%
1718-51-0	Terphenyl-d14	66%		10-157%

(a) Not accredited for this compound at the time of analysis, but all method requirements were followed.



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID:	RA-17 (17-18)	Date Sampled:	04/22/16
Lab Sample ID:	JC19023-1	Date Received:	04/26/16
Matrix:	SO - Soil	Percent Solids:	72.6
Method:	SW846 8015C MOD		
Project:	BMSMC, Building 5 Area, PR		

	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1 ^a	XY064105.D	1	05/03/16	AFL	n/a	n/a	F:GXY2773
Run #2							

	Initial Weight	Final Volume
Run #1	5.34 g	10.0 ml
Run #2		

CAS No.	Compound	Result	RL	MDL	Units	Q
64-17-5	Ethanol	ND	13	2.6	mg/kg	
78-83-1	Isobutyl Alcohol	ND	13	2.6	mg/kg	
67-63-0	Isopropyl Alcohol	ND	13	2.6	mg/kg	
71-23-8	n-Propyl Alcohol	ND	13	2.6	mg/kg	
71-36-3	n-Butyl Alcohol	ND	13	2.6	mg/kg	
67-56-1	Methanol	ND	13	2.6	mg/kg	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
111-27-3	Hexanol	111%		69-121%

(a) Sample was received in a bulk container but was not preserved within 48 hours of sampling. Analysis performed at Accutest Laboratories, Orlando FL.



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID:	RA-17 (17-18)	Date Sampled:	04/22/16
Lab Sample ID:	JC19023-1	Date Received:	04/26/16
Matrix:	SO - Soil	Percent Solids:	72.6
Method:	SW846 8081B SW846 3546		
Project:	BMSMC, Building 5 Area, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	1G122627.D	1	04/29/16	BP	04/28/16	OP93471	G1G3976
Run #2							

Run #	Initial Weight	Final Volume
Run #1	15.4 g	10.0 ml
Run #2		

Pesticide TCL List

CAS No.	Compound	Result	RL	MDL	Units	Q
309-00-2	Aldrin	ND	0.89	0.80	ug/kg	
319-84-6	alpha-BHC	ND	0.89	0.60	ug/kg	
319-85-7	beta-BHC	ND	0.89	0.55	ug/kg	
319-86-8	delta-BHC	ND	0.89	0.35	ug/kg	
58-89-9	gamma-BHC (Lindane)	ND	0.89	0.41	ug/kg	
5103-71-9	alpha-Chlordane	ND	0.89	0.48	ug/kg	
5103-74-2	gamma-Chlordane	ND	0.89	0.68	ug/kg	
60-57-1	Dieldrin	ND	0.89	0.70	ug/kg	
72-54-8	4,4'-DDD	ND	0.89	0.33	ug/kg	
72-55-9	4,4'-DDE	ND	0.89	0.30	ug/kg	
50-29-3	4,4'-DDT	ND	0.89	0.34	ug/kg	
72-20-8	Endrin	ND	0.89	0.32	ug/kg	
1031-07-8	Endosulfan sulfate	ND	0.89	0.51	ug/kg	
7421-93-4	Endrin aldehyde	ND	0.89	0.66	ug/kg	
959-98-8	Endosulfan-I	ND	0.89	0.29	ug/kg	
33213-65-9	Endosulfan-II	ND	0.89	0.85	ug/kg	
76-44-8	Heptachlor	ND	0.89	0.74	ug/kg	
1024-57-3	Heptachlor epoxide	ND	0.89	0.37	ug/kg	
72-43-5	Methoxychlor	ND	1.8	0.50	ug/kg	
53494-70-5	Endrin ketone	ND	0.89	0.47	ug/kg	
8001-35-2	Toxaphene	ND	22	15	ug/kg	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
877-09-8	Tetrachloro-m-xylene	59%		24-136%
877-09-8	Tetrachloro-m-xylene	65%		24-136%
2051-24-3	Decachlorobiphenyl	59%		10-153%
2051-24-3	Decachlorobiphenyl	69%		10-153%



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 3

Client Sample ID:	S-40S (10.5-11.5)	Date Sampled:	04/22/16
Lab Sample ID:	JC19023-2	Date Received:	04/26/16
Matrix:	SO - Soil	Percent Solids:	74.6
Method:	SW846 8270D SW846 3546		
Project:	BMSMC, Building 5 Area, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	5P27908.D	1	05/03/16	AD	04/29/16	OP93478	E5P1416
Run #2							

Run #	Initial Weight	Final Volume
Run #1	30.1 g	1.0 ml
Run #2		

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units	Q
95-57-8	2-Chlorophenol	ND	89	33	ug/kg	
59-50-7	4-Chloro-3-methyl phenol	ND	220	40	ug/kg	
120-83-2	2,4-Dichlorophenol	ND	220	36	ug/kg	
105-67-9	2,4-Dimethylphenol	ND	220	81	ug/kg	
51-28-5	2,4-Dinitrophenol	ND	220	200	ug/kg	
534-52-1	4,6-Dinitro-o-cresol	ND	220	85	ug/kg	
95-48-7	2-Methylphenol	ND	89	65	ug/kg	
	3&4-Methylphenol	ND	89	43	ug/kg	
88-75-5	2-Nitrophenol	ND	220	41	ug/kg	
100-02-7	4-Nitrophenol	ND	450	76	ug/kg	
87-86-5	Pentachlorophenol	ND	220	110	ug/kg	
108-95-2	Phenol	ND	89	33	ug/kg	
58-90-2	2,3,4,6-Tetrachlorophenol	ND	220	42	ug/kg	
95-95-4	2,4,5-Trichlorophenol	ND	220	40	ug/kg	
88-06-2	2,4,6-Trichlorophenol	ND	220	36	ug/kg	
83-32-9	Acenaphthene	ND	45	42	ug/kg	
208-96-8	Acenaphthylene	ND	45	4.7	ug/kg	
98-86-2	Acetophenone	ND	220	7.6	ug/kg	
120-12-7	Anthracene	ND	45	3.8	ug/kg	
1912-24-9	Atrazine	ND	89	18	ug/kg	
56-55-3	Benzo(a)anthracene	ND	45	8.6	ug/kg	
50-32-8	Benzo(a)pyrene	ND	45	9.5	ug/kg	
205-99-2	Benzo(b)fluoranthene	ND	45	9.2	ug/kg	
191-24-2	Benzo(g,h,i)perylene	ND	45	13	ug/kg	
207-08-9	Benzo(k)fluoranthene	ND	45	9.9	ug/kg	
101-55-3	4-Bromophenyl phenyl ether	ND	89	10	ug/kg	
85-68-7	Butyl benzyl phthalate	ND	89	24	ug/kg	
92-52-4	1,1'-Biphenyl	ND	89	8.2	ug/kg	
100-52-7	Benzaldehyde	ND	220	11	ug/kg	
91-58-7	2-Chloronaphthalene	ND	89	6.4	ug/kg	
106-47-8	4-Chloroaniline	ND	220	12	ug/kg	
86-74-8	Carbazole	ND	89	4.9	ug/kg	



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

Report of Analysis

Client Sample ID: S-40S (10.5-11.5)
 Lab Sample ID: JC19023-2
 Matrix: SO - Soil
 Method: SW846 8270D SW846 3546
 Project: BSMC, Building 5 Area, PR

Date Sampled: 04/22/16
 Date Received: 04/26/16
 Percent Solids: 74.6

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units	Q
105-60-2	Caprolactam	ND	89	29	ug/kg	
218-01-9	Chrysene	ND	45	7.2	ug/kg	
111-91-1	bis(2-Chloroethoxy)methane	ND	89	10	ug/kg	
111-44-4	bis(2-Chloroethyl)ether	ND	89	19	ug/kg	
108-60-1	bis(2-Chloroisopropyl)ether	ND	89	10	ug/kg	
7005-72-3	4-Chlorophenyl phenyl ether	ND	89	8.4	ug/kg	
121-14-2	2,4-Dinitrotoluene	ND	45	8.4	ug/kg	
606-20-2	2,6-Dinitrotoluene	ND	45	11	ug/kg	
91-94-1	3,3'-Dichlorobenzidine	ND	89	29	ug/kg	
53-70-3	Dibenzo(a,h)anthracene	ND	45	16	ug/kg	
132-64-9	Dibenzofuran	ND	89	6.2	ug/kg	
84-74-2	Di-n-butyl phthalate	ND	89	5.3	ug/kg	
117-84-0	Di-n-octyl phthalate	ND	89	6.0	ug/kg	
84-66-2	Diethyl phthalate	ND	89	5.7	ug/kg	
131-11-3	Dimethyl phthalate	ND	89	6.4	ug/kg	
117-81-7	bis(2-Ethylhexyl)phthalate	ND	89	16	ug/kg	
206-44-0	Fluoranthene	ND	45	5.4	ug/kg	
86-73-7	Fluorene	ND	45	5.3	ug/kg	
118-74-1	Hexachlorobenzene	ND	89	8.8	ug/kg	
87-68-3	Hexachlorobutadiene	ND	45	12	ug/kg	
77-47-4	Hexachlorocyclopentadiene	ND	450	71	ug/kg	
67-72-1	Hexachloroethane	ND	220	14	ug/kg	
193-39-5	Indeno(1,2,3-cd)pyrene	ND	45	23	ug/kg	
78-59-1	Isophorone	ND	89	8.3	ug/kg	
90-12-0	1-Methylnaphthalene	ND	89	7.3	ug/kg	
91-57-6	2-Methylnaphthalene	ND	89	8.3	ug/kg	
88-74-4	2-Nitroaniline	ND	220	10	ug/kg	
99-09-2	3-Nitroaniline	ND	220	13	ug/kg	
100-01-6	4-Nitroaniline	ND	220	15	ug/kg	
98-95-3	Nitrobenzene	ND	89	14	ug/kg	
621-64-7	N-Nitroso-di-n-propylamine	ND	89	13	ug/kg	
86-30-6	N-Nitrosodiphenylamine	ND	220	23	ug/kg	
85-01-8	Phenanthrene	ND	45	4.9	ug/kg	
129-00-0	Pyrene	ND	45	5.6	ug/kg	
95-94-3	1,2,4,5-Tetrachlorobenzene	ND	220	11	ug/kg	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
367-12-4	2-Fluorophenol	72%		30-106%
4165-62-2	Phenol-d5	72%		30-106%

ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound



Report of Analysis

Client Sample ID:	S-40S (10.5-11.5)	Date Sampled:	04/22/16
Lab Sample ID:	JC19023-2	Date Received:	04/26/16
Matrix:	SO - Soil	Percent Solids:	74.6
Method:	SW846 8270D SW846 3546		
Project:	BMSMC, Building 5 Area, PR		

ABN TCL Special List

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
118-79-6	2,4,6-Tribromophenol	76%		24-140%
4165-60-0	Nitrobenzene-d5	77%		26-122%
321-60-8	2-Fluorobiphenyl	80%		36-112%
1718-51-0	Terphenyl-d14	90%		36-132%



ND = Not detected MDL = Method Detection Limit
RL = Reporting Limit
E = Indicates value exceeds calibration range

J = Indicates an estimated value
B = Indicates analyte found in associated method blank
N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID:	S-40S (10.5-11.5)	Date Sampled:	04/22/16
Lab Sample ID:	JC19023-2	Date Received:	04/26/16
Matrix:	SO - Soil	Percent Solids:	74.6
Method:	SW846 8270D BY SIM SW846 3546		
Project:	BMSMC, Building 5 Area, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	4M65190.D	1	05/09/16	LK	04/29/16	OP93478A	E4M2910
Run #2							

Run #	Initial Weight	Final Volume
Run #1	30.1 g	1.0 ml
Run #2		

CAS No.	Compound	Result	RL	MDL	Units	Q
123-91-1	1,4-Dioxane ^a	ND	4.5	0.90	ug/kg	
91-20-3	Naphthalene	ND	4.5	0.54	ug/kg	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
4165-60-0	Nitrobenzene-d5	89%		15-138%
321-60-8	2-Fluorobiphenyl	77%		12-148%
1718-51-0	Terphenyl-d14	86%		10-157%

(a) Not accredited for this compound at the time of analysis, but all method requirements were followed.



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID:	S-40S (10.5-11.5)	Date Sampled:	04/22/16
Lab Sample ID:	JC19023-2	Date Received:	04/26/16
Matrix:	SO - Soil	Percent Solids:	74.6
Method:	SW846 8015C MOD		
Project:	BMSMC, Building 5 Area, PR		

	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1 ^a	XY064106.D	1	05/03/16	AFL	n/a	n/a	F:GXY2773
Run #2							

	Initial Weight	Final Volume
Run #1	4.96 g	10.0 ml
Run #2		

CAS No.	Compound	Result	RL	MDL	Units	Q
64-17-5	Ethanol	ND	14	2.7	mg/kg	
78-83-1	Isobutyl Alcohol	ND	14	2.7	mg/kg	
67-63-0	Isopropyl Alcohol	ND	14	2.7	mg/kg	
71-23-8	n-Propyl Alcohol	ND	14	2.7	mg/kg	
71-36-3	n-Butyl Alcohol	ND	14	2.7	mg/kg	
67-56-1	Methanol	ND	14	2.7	mg/kg	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
111-27-3	Hexanol	137% ^b		69-121%

(a) Sample was received in a bulk container but was not preserved within 48 hours of sampling. Analysis performed at Accutest Laboratories, Orlando FL.

(b) Outside control limits; however, sample is ND.



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID:	S-40S (10.5-11.5)	Date Sampled:	04/22/16
Lab Sample ID:	JC19023-2	Date Received:	04/26/16
Matrix:	SO - Soil	Percent Solids:	74.6
Method:	SW846 8081B SW846 3546		
Project:	BMSMC, Building 5 Area, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	1G122628.D	1	04/29/16	BP	04/28/16	OP93471	G1G3976
Run #2							

Run #	Initial Weight	Final Volume
Run #1	16.3 g	10.0 ml
Run #2		

Pesticide TCL List

CAS No.	Compound	Result	RL	MDL	Units	Q
309-00-2	Aldrin	ND	0.82	0.73	ug/kg	
319-84-6	alpha-BHC	ND	0.82	0.55	ug/kg	
319-85-7	beta-BHC	ND	0.82	0.51	ug/kg	
319-86-8	delta-BHC	ND	0.82	0.32	ug/kg	
58-89-9	gamma-BHC (Lindane)	ND	0.82	0.37	ug/kg	
5103-71-9	alpha-Chlordane	ND	0.82	0.44	ug/kg	
5103-74-2	gamma-Chlordane	ND	0.82	0.63	ug/kg	
60-57-1	Dieldrin	ND	0.82	0.64	ug/kg	
72-54-8	4,4'-DDD	ND	0.82	0.30	ug/kg	
72-55-9	4,4'-DDE	ND	0.82	0.27	ug/kg	
50-29-3	4,4'-DDT	ND	0.82	0.31	ug/kg	
72-20-8	Endrin	ND	0.82	0.29	ug/kg	
1031-07-8	Endosulfan sulfate	ND	0.82	0.47	ug/kg	
7421-93-4	Endrin aldehyde	ND	0.82	0.61	ug/kg	
959-98-8	Endosulfan-I	ND	0.82	0.27	ug/kg	
33213-65-9	Endosulfan-II	ND	0.82	0.78	ug/kg	
76-44-8	Heptachlor	ND	0.82	0.68	ug/kg	
1024-57-3	Heptachlor epoxide	ND	0.82	0.34	ug/kg	
72-43-5	Methoxychlor	ND	1.6	0.46	ug/kg	
53494-70-5	Endrin ketone	ND	0.82	0.43	ug/kg	
8001-35-2	Toxaphene	ND	21	14	ug/kg	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
877-09-8	Tetrachloro-m-xylene	80%		24-136%
877-09-8	Tetrachloro-m-xylene	84%		24-136%
2051-24-3	Decachlorobiphenyl	76%		10-153%
2051-24-3	Decachlorobiphenyl	82%		10-153%



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 3

Client Sample ID:	BPEB-6	Date Sampled:	04/22/16
Lab Sample ID:	JC19023-3	Date Received:	04/26/16
Matrix:	AQ - Equipment Blank	Percent Solids:	n/a
Method:	SW846 8270D SW846 3510C		
Project:	BMSMC, Building 5 Area, PR		

	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	2P58712.D	1	04/29/16	RL	04/28/16	OP93453	E2P2566
Run #2							

	Initial Volume	Final Volume
Run #1	1000 ml	1.0 ml
Run #2		

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units	Q
95-57-8	2-Chlorophenol	ND	5.0	0.82	ug/l	
59-50-7	4-Chloro-3-methyl phenol	ND	5.0	0.89	ug/l	
120-83-2	2,4-Dichlorophenol	ND	2.0	1.3	ug/l	
105-67-9	2,4-Dimethylphenol	ND	5.0	2.4	ug/l	
51-28-5	2,4-Dinitrophenol	ND	10	1.6	ug/l	
534-52-1	4,6-Dinitro-o-cresol	ND	5.0	1.3	ug/l	
95-48-7	2-Methylphenol	ND	2.0	0.89	ug/l	
	3&4-Methylphenol	ND	2.0	0.88	ug/l	
88-75-5	2-Nitrophenol	ND	5.0	0.96	ug/l	
100-02-7	4-Nitrophenol	ND	10	1.2	ug/l	
87-86-5	Pentachlorophenol	ND	5.0	1.4	ug/l	
108-95-2	Phenol	ND	2.0	0.39	ug/l	
58-90-2	2,3,4,6-Tetrachlorophenol	ND	5.0	1.5	ug/l	
95-95-4	2,4,5-Trichlorophenol	ND	5.0	1.3	ug/l	
88-06-2	2,4,6-Trichlorophenol	ND	5.0	0.92	ug/l	
83-32-9	Acenaphthene	ND	1.0	0.19	ug/l	
208-96-8	Acenaphthylene	ND	1.0	0.14	ug/l	
98-86-2	Acetophenone	ND	2.0	0.21	ug/l	
120-12-7	Anthracene	ND	1.0	0.21	ug/l	
1912-24-9	Atrazine	ND	2.0	0.45	ug/l	
100-52-7	Benzaldehyde	ND	5.0	0.29	ug/l	
56-55-3	Benzo(a)anthracene	ND	1.0	0.20	ug/l	
50-32-8	Benzo(a)pyrene	ND	1.0	0.21	ug/l	
205-99-2	Benzo(b)fluoranthene	ND	1.0	0.21	ug/l	
191-24-2	Benzo(g,h,i)perylene	ND	1.0	0.34	ug/l	
207-08-9	Benzo(k)fluoranthene	ND	1.0	0.21	ug/l	
101-55-3	4-Bromophenyl phenyl ether	ND	2.0	0.40	ug/l	
85-68-7	Butyl benzyl phthalate	ND	2.0	0.46	ug/l	
92-52-4	1,1'-Biphenyl	ND	1.0	0.21	ug/l	
91-58-7	2-Chloronaphthalene	ND	2.0	0.24	ug/l	
106-47-8	4-Chloroaniline	ND	5.0	0.34	ug/l	
86-74-8	Carbazole	ND	1.0	0.23	ug/l	

ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound



Report of Analysis

Client Sample ID: BPEB-6
 Lab Sample ID: JC19023-3
 Matrix: AQ - Equipment Blank
 Method: SW846 8270D SW846 3510C
 Project: BMSMC, Building 5 Area, PR

Date Sampled: 04/22/16
 Date Received: 04/26/16
 Percent Solids: n/a

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units	Q
105-60-2	Caprolactam	ND	2.0	0.65	ug/l	
218-01-9	Chrysene	ND	1.0	0.18	ug/l	
111-91-1	bis(2-Chloroethoxy)methane	ND	2.0	0.28	ug/l	
111-44-4	bis(2-Chloroethyl)ether	ND	2.0	0.25	ug/l	
108-60-1	bis(2-Chloroisopropyl)ether	ND	2.0	0.40	ug/l	
7005-72-3	4-Chlorophenyl phenyl ether	ND	2.0	0.37	ug/l	
121-14-2	2,4-Dinitrotoluene	ND	1.0	0.55	ug/l	
606-20-2	2,6-Dinitrotoluene	ND	1.0	0.48	ug/l	
91-94-1	3,3'-Dichlorobenzidine	ND	2.0	0.51	ug/l	
53-70-3	Dibenzo(a,h)anthracene	ND	1.0	0.33	ug/l	
132-64-9	Dibenzofuran	ND	5.0	0.22	ug/l	
84-74-2	Di-n-butyl phthalate	1.8	2.0	0.50	ug/l	JB
117-84-0	Di-n-octyl phthalate	ND	2.0	0.23	ug/l	
84-66-2	Diethyl phthalate	ND	2.0	0.26	ug/l	
131-11-3	Dimethyl phthalate	ND	2.0	0.22	ug/l	
117-81-7	bis(2-Ethylhexyl)phthalate	ND	2.0	1.7	ug/l	
206-44-0	Fluoranthene	ND	1.0	0.17	ug/l	
86-73-7	Fluorene	ND	1.0	0.17	ug/l	
118-74-1	Hexachlorobenzene	ND	1.0	0.33	ug/l	
87-68-3	Hexachlorobutadiene	ND	1.0	0.49	ug/l	
77-47-4	Hexachlorocyclopentadiene	ND	10	2.8	ug/l	
67-72-1	Hexachloroethane	ND	2.0	0.39	ug/l	
193-39-5	Indeno(1,2,3-cd)pyrene	ND	1.0	0.33	ug/l	
78-59-1	Isophorone	ND	2.0	0.28	ug/l	
90-12-0	1-Methylnaphthalene	ND	1.0	0.26	ug/l	
91-57-6	2-Methylnaphthalene	ND	1.0	0.21	ug/l	
88-74-4	2-Nitroaniline	ND	5.0	0.28	ug/l	
99-09-2	3-Nitroaniline	ND	5.0	0.39	ug/l	
100-01-6	4-Nitroaniline	ND	5.0	0.44	ug/l	
98-95-3	Nitrobenzene	ND	2.0	0.64	ug/l	
621-64-7	N-Nitroso-di-n-propylamine	ND	2.0	0.48	ug/l	
86-30-6	N-Nitrosodiphenylamine	ND	5.0	0.22	ug/l	
85-01-8	Phenanthrene	ND	1.0	0.18	ug/l	
129-00-0	Pyrene	ND	1.0	0.22	ug/l	
95-94-3	1,2,4,5-Tetrachlorobenzene	ND	2.0	0.37	ug/l	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
367-12-4	2-Fluorophenol	33%		14-88%
4165-62-2	Phenol-d5	23%		10-110%

ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound



Report of Analysis

Client Sample ID:	BPEB-6	Date Sampled:	04/22/16
Lab Sample ID:	JC19023-3	Date Received:	04/26/16
Matrix:	AQ - Equipment Blank	Percent Solids:	n/a
Method:	SW846 8270D SW846 3510C		
Project:	BMSMC, Building 5 Area, PR		

ABN TCL Special List

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
118-79-6	2,4,6-Tribromophenol	79%		39-149%
4165-60-0	Nitrobenzene-d5	63%		32-128%
321-60-8	2-Fluorobiphenyl	64%		35-119%
1718-51-0	Terphenyl-d14	86%		10-126%



ND = Not detected MDL = Method Detection Limit
RL = Reporting Limit
E = Indicates value exceeds calibration range

J = Indicates an estimated value
B = Indicates analyte found in associated method blank
N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID:	BPEB-6	Date Sampled:	04/22/16
Lab Sample ID:	JC19023-3	Date Received:	04/26/16
Matrix:	AQ - Equipment Blank	Percent Solids:	n/a
Method:	SW846 8270D BY SIM SW846 3510C		
Project:	BMSMC, Building 5 Area, PR		

	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	3M60957.D	1	04/29/16	JJ	04/28/16	OP93453A	E3M2866
Run #2							

	Initial Volume	Final Volume
Run #1	1000 ml	1.0 ml
Run #2		

CAS No.	Compound	Result	RL	MDL	Units	Q
91-20-3	Naphthalene	ND	0.10	0.029	ug/l	
123-91-1	1,4-Dioxane	ND	0.10	0.049	ug/l	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
4165-60-0	Nitrobenzene-d5	57%		24-125%
321-60-8	2-Fluorobiphenyl	69%		19-127%
1718-51-0	Terphenyl-d14	92%		10-119%



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID: BPEB-6		Date Sampled: 04/22/16
Lab Sample ID: JC19023-3		Date Received: 04/26/16
Matrix: AQ - Equipment Blank		Percent Solids: n/a
Method: SW846 8015C		
Project: BMSMC, Building 5 Area, PR		

Run #1 ^a	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #2	XY064092.D	1	05/02/16	AFL	n/a	n/a	F:GXY2772

CAS No.	Compound	Result	RL	MDL	Units	Q
64-17-5	Ethanol	ND	5.0	1.0	mg/l	
78-83-1	Isobutyl Alcohol	ND	5.0	1.0	mg/l	
67-63-0	Isopropyl Alcohol	ND	5.0	1.0	mg/l	
71-23-8	n-Propyl Alcohol	ND	5.0	1.0	mg/l	
71-36-3	n-Butyl Alcohol	ND	5.0	1.0	mg/l	
67-56-1	Methanol	ND	5.0	1.0	mg/l	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
111-27-3	Hexanol	91%		73-123%

(a) Analysis performed at Accutest Laboratories, Orlando FL.



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID:	BPEB-6	Date Sampled:	04/22/16
Lab Sample ID:	JC19023-3	Date Received:	04/26/16
Matrix:	AQ - Equipment Blank	Percent Solids:	n/a
Method:	SW846 8081B SW846 3510C		
Project:	BMSMC, Building 5 Area, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	6G34452.D	1	04/27/16	BP	04/26/16	OP93399	G6G997
Run #2							

Run #	Initial Volume	Final Volume
Run #1	300 ml	2.0 ml
Run #2		

Pesticide TCL List

CAS No.	Compound	Result	RL	MDL	Units	Q
309-00-2	Aldrin	ND	0.0067	0.0040	ug/l	
319-84-6	alpha-BHC	ND	0.0067	0.0040	ug/l	
319-85-7	beta-BHC	ND	0.0067	0.0038	ug/l	
319-86-8	delta-BHC	ND	0.0067	0.0030	ug/l	
58-89-9	gamma-BHC (Lindane)	ND	0.0067	0.0019	ug/l	
5103-71-9	alpha-Chlordane	ND	0.0067	0.0031	ug/l	
5103-74-2	gamma-Chlordane	ND	0.0067	0.0031	ug/l	
60-57-1	Dieldrin	ND	0.0067	0.0024	ug/l	
72-54-8	4,4'-DDD	ND	0.0067	0.0025	ug/l	
72-55-9	4,4'-DDE	ND	0.0067	0.0041	ug/l	
50-29-3	4,4'-DDT	ND	0.0067	0.0033	ug/l	
72-20-8	Endrin	ND	0.0067	0.0034	ug/l	
1031-07-8	Endosulfan sulfate	ND	0.0067	0.0035	ug/l	
7421-93-4	Endrin aldehyde	ND	0.0067	0.0034	ug/l	
53494-70-5	Endrin ketone	ND	0.0067	0.0034	ug/l	
959-98-8	Endosulfan-I	ND	0.0067	0.0033	ug/l	
33213-65-9	Endosulfan-II	ND	0.0067	0.0029	ug/l	
76-44-8	Heptachlor	ND	0.0067	0.0025	ug/l	
1024-57-3	Heptachlor epoxide	ND	0.0067	0.0044	ug/l	
72-43-5	Methoxychlor	ND	0.013	0.0038	ug/l	
8001-35-2	Toxaphene	ND	0.17	0.12	ug/l	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
877-09-8	Tetrachloro-m-xylene	96%		26-132%
877-09-8	Tetrachloro-m-xylene	97%		26-132%
2051-24-3	Decachlorobiphenyl	61%		10-118%
2051-24-3	Decachlorobiphenyl	67%		10-118%



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 3

Client Sample ID:	S-39 (11-12)	Date Sampled:	04/25/16
Lab Sample ID:	JC19023-4	Date Received:	04/26/16
Matrix:	SO - Soil	Percent Solids:	75.4
Method:	SW846 8270D SW846 3546		
Project:	BMSMC, Building 5 Area, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	5P27909.D	1	05/03/16	AD	04/29/16	OP93478	E5P1416
Run #2							

Run #	Initial Weight	Final Volume
Run #1	31.3 g	1.0 ml
Run #2		

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units	Q
95-57-8	2-Chlorophenol	ND	85	31	ug/kg	
59-50-7	4-Chloro-3-methyl phenol	ND	210	38	ug/kg	
120-83-2	2,4-Dichlorophenol	ND	210	34	ug/kg	
105-67-9	2,4-Dimethylphenol	ND	210	78	ug/kg	
51-28-5	2,4-Dinitrophenol	ND	210	190	ug/kg	
534-52-1	4,6-Dinitro-o-cresol	ND	210	81	ug/kg	
95-48-7	2-Methylphenol	ND	85	61	ug/kg	
	3&4-Methylphenol	ND	85	41	ug/kg	
88-75-5	2-Nitrophenol	ND	210	39	ug/kg	
100-02-7	4-Nitrophenol	ND	420	72	ug/kg	
87-86-5	Pentachlorophenol	ND	210	100	ug/kg	
108-95-2	Phenol	ND	85	32	ug/kg	
58-90-2	2,3,4,6-Tetrachlorophenol	ND	210	40	ug/kg	
95-95-4	2,4,5-Trichlorophenol	ND	210	38	ug/kg	
88-06-2	2,4,6-Trichlorophenol	ND	210	34	ug/kg	
83-32-9	Acenaphthene	ND	42	40	ug/kg	
208-96-8	Acenaphthylene	ND	42	4.4	ug/kg	
98-86-2	Acetophenone	ND	210	7.2	ug/kg	
120-12-7	Anthracene	72.1	42	3.6	ug/kg	
1912-24-9	Atrazine	ND	85	17	ug/kg	
56-55-3	Benzo(a)anthracene	ND	42	8.2	ug/kg	
50-32-8	Benzo(a)pyrene	ND	42	9.0	ug/kg	
205-99-2	Benzo(b)fluoranthene	ND	42	8.7	ug/kg	
191-24-2	Benzo(g,h,i)perylene	ND	42	13	ug/kg	
207-08-9	Benzo(k)fluoranthene	ND	42	9.4	ug/kg	
101-55-3	4-Bromophenyl phenyl ether	ND	85	9.7	ug/kg	
85-68-7	Butyl benzyl phthalate	ND	85	23	ug/kg	
92-52-4	1,1'-Biphenyl	ND	85	7.8	ug/kg	
100-52-7	Benzaldehyde	ND	210	11	ug/kg	
91-58-7	2-Chloronaphthalene	ND	85	6.1	ug/kg	
106-47-8	4-Chloroaniline	ND	210	11	ug/kg	
86-74-8	Carbazole	ND	85	4.7	ug/kg	

ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound



Report of Analysis

Client Sample ID:	S-39 (11-12)	Date Sampled:	04/25/16
Lab Sample ID:	JC19023-4	Date Received:	04/26/16
Matrix:	SO - Soil	Percent Solids:	75.4
Method:	SW846 8270D SW846 3546		
Project:	BMSMC, Building 5 Area, PR		

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units	Q
105-60-2	Caprolactam	ND	85	27	ug/kg	
218-01-9	Chrysene	ND	42	6.8	ug/kg	
111-91-1	bis(2-Chloroethoxy)methane	ND	85	9.6	ug/kg	
111-44-4	bis(2-Chloroethyl)ether	ND	85	18	ug/kg	
108-60-1	bis(2-Chloroisopropyl)ether	ND	85	9.7	ug/kg	
7005-72-3	4-Chlorophenyl phenyl ether	ND	85	8.0	ug/kg	
121-14-2	2,4-Dinitrotoluene	ND	42	8.0	ug/kg	
606-20-2	2,6-Dinitrotoluene	ND	42	11	ug/kg	
91-94-1	3,3'-Dichlorobenzidine	ND	85	28	ug/kg	
53-70-3	Dibenzo(a,h)anthracene	ND	42	15	ug/kg	
132-64-9	Dibenzofuran	ND	85	5.9	ug/kg	
84-74-2	Di-n-butyl phthalate	ND	85	5.0	ug/kg	
117-84-0	Di-n-octyl phthalate	ND	85	5.7	ug/kg	
84-66-2	Diethyl phthalate	ND	85	5.4	ug/kg	
131-11-3	Dimethyl phthalate	ND	85	6.1	ug/kg	
117-81-7	bis(2-Ethylhexyl)phthalate	ND	85	15	ug/kg	
206-44-0	Fluoranthene	ND	42	5.2	ug/kg	
86-73-7	Fluorene	ND	42	5.0	ug/kg	
118-74-1	Hexachlorobenzene	ND	85	8.3	ug/kg	
87-68-3	Hexachlorobutadiene	ND	42	11	ug/kg	
77-47-4	Hexachlorocyclopentadiene	ND	420	67	ug/kg	
67-72-1	Hexachloroethane	ND	210	14	ug/kg	
193-39-5	Indeno(1,2,3-cd)pyrene	ND	42	22	ug/kg	
78-59-1	Isophorone	ND	85	7.9	ug/kg	
90-12-0	1-Methylnaphthalene	ND	85	6.9	ug/kg	
91-57-6	2-Methylnaphthalene	ND	85	7.9	ug/kg	
88-74-4	2-Nitroaniline	ND	210	9.6	ug/kg	
99-09-2	3-Nitroaniline	ND	210	12	ug/kg	
100-01-6	4-Nitroaniline	ND	210	14	ug/kg	
98-95-3	Nitrobenzene	ND	85	13	ug/kg	
621-64-7	N-Nitroso-di-n-propylamine	ND	85	13	ug/kg	
86-30-6	N-Nitrosodiphenylamine	ND	210	22	ug/kg	
85-01-8	Phenanthrene	ND	42	4.7	ug/kg	
129-00-0	Pyrene	ND	42	5.3	ug/kg	
95-94-3	1,2,4,5-Tetrachlorobenzene	ND	210	10	ug/kg	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
367-12-4	2-Fluorophenol	80%		30-106%
4165-62-2	Phenol-d5	80%		30-106%

ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound



Report of Analysis

Client Sample ID: S-39 (11-12)	
Lab Sample ID: JC19023-4	Date Sampled: 04/25/16
Matrix: SO - Soil	Date Received: 04/26/16
Method: SW846 8270D SW846 3546	Percent Solids: 75.4
Project: BMSMC, Building 5 Area, PR	

ABN TCL Special List

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
118-79-6	2,4,6-Tribromophenol	72%		24-140%
4165-60-0	Nitrobenzene-d5	32%		26-122%
321-60-8	2-Fluorobiphenyl	80%		36-112%
1718-51-0	Terphenyl-d14	89%		36-132%



ND = Not detected MDL = Method Detection Limit
RL = Reporting Limit
E = Indicates value exceeds calibration range

J = Indicates an estimated value
B = Indicates analyte found in associated method blank
N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID:	S-39 (11-12)	Date Sampled:	04/25/16
Lab Sample ID:	JC19023-4	Date Received:	04/26/16
Matrix:	SO - Soil	Percent Solids:	75.4
Method:	SW846 8270D BY SIM SW846 3546		
Project:	BMSMC, Building 5 Area, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	4M65191.D	1	05/09/16	LK	04/29/16	OP93478A	E4M2910
Run #2							

Run #	Initial Weight	Final Volume
Run #1	31.3 g	1.0 ml
Run #2		

CAS No.	Compound	Result	RL	MDL	Units	Q
123-91-1	1,4-Dioxane ^a	ND	4.2	0.85	ug/kg	
91-20-3	Naphthalene	ND	4.2	0.52	ug/kg	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
4165-60-0	Nitrobenzene-d5	35%		15-138%
321-60-8	2-Fluorobiphenyl	79%		12-148%
1718-51-0	Terphenyl-d14	88%		10-157%

(a) Not accredited for this compound at the time of analysis, but all method requirements were followed.



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID:	S-39 (11-12)	Date Sampled:	04/25/16
Lab Sample ID:	JC19023-4	Date Received:	04/26/16
Matrix:	SO - Soil	Percent Solids:	75.4
Method:	SW846 8015C MOD		
Project:	BMSMC, Building 5 Area, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1 ^a	XY064107.D	1	05/03/16	AFL	n/a	n/a	F:GXY2773
Run #2							

Run #	Initial Weight	Final Volume
Run #1	5.12 g	10.0 ml
Run #2		

CAS No.	Compound	Result	RL	MDL	Units	Q
64-17-5	Ethanol	ND	13	2.6	mg/kg	
78-83-1	Isobutyl Alcohol	ND	13	2.6	mg/kg	
67-63-0	Isopropyl Alcohol	ND	13	2.6	mg/kg	
71-23-8	n-Propyl Alcohol	ND	13	2.6	mg/kg	
71-36-3	n-Butyl Alcohol	ND	13	2.6	mg/kg	
67-56-1	Methanol	ND	13	2.6	mg/kg	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
111-27-3	Hexanol	96%		69-121%

(a) Sample was received in a bulk container but was not preserved within 48 hours of sampling. Analysis performed at Accutest Laboratories, Orlando FL.



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID:	S-39 (11-12)	Date Sampled:	04/25/16
Lab Sample ID:	JC19023-4	Date Received:	04/26/16
Matrix:	SO - Soil	Percent Solids:	75.4
Method:	SW846 8081B SW846 3546		
Project:	BMSMC, Building 5 Area, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	1G122629.D	1	04/29/16	BP	04/28/16	OP93471	G1G3976
Run #2							

Run #	Initial Weight	Final Volume
Run #1	15.7 g	10.0 ml
Run #2		

Pesticide TCL List

CAS No.	Compound	Result	RL	MDL	Units	Q
309-00-2	Aldrin	ND	0.84	0.75	ug/kg	
319-84-6	alpha-BHC	ND	0.84	0.56	ug/kg	
319-85-7	beta-BHC	ND	0.84	0.52	ug/kg	
319-86-8	delta-BHC	ND	0.84	0.33	ug/kg	
58-89-9	gamma-BHC (Lindane)	ND	0.84	0.38	ug/kg	
5103-71-9	alpha-Chlordane	ND	0.84	0.45	ug/kg	
5103-74-2	gamma-Chlordane	ND	0.84	0.64	ug/kg	
60-57-1	Dieldrin	ND	0.84	0.66	ug/kg	
72-54-8	4,4'-DDD	ND	0.84	0.31	ug/kg	
72-55-9	4,4'-DDE	ND	0.84	0.28	ug/kg	
50-29-3	4,4'-DDT	ND	0.84	0.32	ug/kg	
72-20-8	Endrin	ND	0.84	0.30	ug/kg	
1031-07-8	Endosulfan sulfate	ND	0.84	0.48	ug/kg	
7421-93-4	Endrin aldehyde	ND	0.84	0.63	ug/kg	
959-98-8	Endosulfan-I	ND	0.84	0.28	ug/kg	
33213-65-9	Endosulfan-II	ND	0.84	0.80	ug/kg	
76-44-8	Heptachlor	ND	0.84	0.69	ug/kg	
1024-57-3	Heptachlor epoxide	ND	0.84	0.35	ug/kg	
72-43-5	Methoxychlor	ND	1.7	0.47	ug/kg	
53494-70-5	Endrin ketone	ND	0.84	0.44	ug/kg	
8001-35-2	Toxaphene	ND	21	15	ug/kg	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
877-09-8	Tetrachloro-m-xylene	61%		24-136%
877-09-8	Tetrachloro-m-xylene	58%		24-136%
2051-24-3	Decachlorobiphenyl	58%		10-153%
2051-24-3	Decachlorobiphenyl	57%		10-153%



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SD - Ereignis
SEB

CHAIN OF CUSTODY

PAGE 1 OF 1

ACCUTEST - N₁

SGS Accountant - Dayton
2225 Route 130, Dayton, NJ 05410
TEL 732-329-0200 FAX 732-329-3459-3480
www.accountant.com

801.219536143

Terms Order Control 6

5 JUL 1962

[illegible]

JC19023: Chain of Custody
Page 1 of 2

EXECUTIVE NARRATIVE

SDG No: **JC19023** Laboratory: **Accutest, New Jersey**
Analysis: **SW846-8081B** Number of Samples: **4**

Location: **BMSMC, Building 5 Area**
Humacao, PR

SUMMARY: Three (3) soil samples and one (1) equipment blank sample were analyzed for selected pesticides following method SW846-8081B. The sample results were assessed according to USEPA data validation guidance documents in the following order of precedence *Hazardous Waste Support Section SOP No. HW-36A, Revision 0, June, 2015. SOM02.2. Pesticide Data Validation*. The QC criteria and data validation actions listed on the data review worksheets are from the primary guidance document, unless otherwise noted.

Results are valid and can be used for decision making purposes.

Critical issues: **None**
Major: **None**
Minor: **None**

Critical findings: **None**
Major findings: **None**
Minor findings: 1. Blank spike/blank spike duplicate RPD outside laboratory control limits. No action taken, professional judgment.

COMMENTS: Results are valid and can be used for decision making purposes.

Reviewers Name: **Rafael Infante**
Chemist License 1888

Signature:



Date: **May 17, 2016**

SAMPLE ORGANIC DATA SAMPLE SUMMARY

Sample ID: JC19023-1

Sample location: BMSMC Building 5 Area

Sampling date: 22-Apr-16

Matrix: Soil

METHOD: 8081B

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Aldrin	0.89	ug/kg	1	-	U	Yes
alpha-BHC	0.89	ug/kg	1	-	U	Yes
beta-BHC	0.89	ug/kg	1	-	U	Yes
delta-BHC	0.89	ug/kg	1	-	U	Yes
gamma-BHC (Lindane)	0.89	ug/kg	1	-	U	Yes
alpha-Chlordane	0.89	ug/kg	1	-	U	Yes
gamma-Chlordane	0.89	ug/kg	1	-	U	Yes
Dieldrin	0.89	ug/kg	1	-	U	Yes
4,4'-DDD	0.89	ug/kg	1	-	U	Yes
4,4'-DDE	0.89	ug/kg	1	-	U	Yes
4,4'-DDT	0.89	ug/kg	1	-	U	Yes
Endrin	0.89	ug/kg	1	-	U	Yes
Endosulfan sulfate	0.89	ug/kg	1	-	U	Yes
Endrin aldehyde	0.89	ug/kg	1	-	U	Yes
Endosulfan-I	0.89	ug/kg	1	-	U	Yes
Endosulfan-II	0.89	ug/kg	1	-	U	Yes
Heptachlor	0.89	ug/kg	1	-	U	Yes
Heptachlor epoxide	0.89	ug/kg	1	-	U	Yes
Methoxychlor	1.8	ug/kg	1	-	U	Yes
Endrin ketone	0.89	ug/kg	1	-	U	Yes
Toxaphene	22	ug/kg	1	-	U	Yes

Sample ID: JC19023-2
Sample location: BMSMC Building 5 Area
Sampling date: 22-Apr-16
Matrix: Soil

METHOD: 8081B

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Aldrin	0.82	ug/kg	1	-	U	Yes
alpha-BHC	0.82	ug/kg	1	-	U	Yes
beta-BHC	0.82	ug/kg	1	-	U	Yes
delta-BHC	0.82	ug/kg	1	-	U	Yes
gamma-BHC (Lindane)	0.82	ug/kg	1	-	U	Yes
alpha-Chlordane	0.82	ug/kg	1	-	U	Yes
gamma-Chlordane	0.82	ug/kg	1	-	U	Yes
Dieldrin	0.82	ug/kg	1	-	U	Yes
4,4'-DDD	0.82	ug/kg	1	-	U	Yes
4,4'-DDE	0.82	ug/kg	1	-	U	Yes
4,4'-DDT	0.82	ug/kg	1	-	U	Yes
Endrin	0.82	ug/kg	1	-	U	Yes
Endosulfan sulfate	0.82	ug/kg	1	-	U	Yes
Endrin aldehyde	0.82	ug/kg	1	-	U	Yes
Endosulfan-I	0.82	ug/kg	1	-	U	Yes
Endosulfan-II	0.82	ug/kg	1	-	U	Yes
Heptachlor	0.82	ug/kg	1	-	U	Yes
Heptachlor epoxide	0.82	ug/kg	1	-	U	Yes
Methoxychlor	1.6	ug/kg	1	-	U	Yes
Endrin ketone	0.82	ug/kg	1	-	U	Yes
Toxaphene	21	ug/kg	1	-	U	Yes

Sample ID: JC19214-3
Sample location: BMSMC Building 5 Area
Sampling date: 22-Apr-16
Matrix: AQ - Equipment Blank

METHOD: 8081B

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Aldrin	0.0067	ug/L	1	-	U	Yes
alpha-BHC	0.0067	ug/L	1	-	U	Yes
beta-BHC	0.0067	ug/L	1	-	U	Yes
delta-BHC	0.0067	ug/L	1	-	U	Yes
gamma-BHC (Lindane)	0.0067	ug/L	1	-	U	Yes
alpha-Chlordane	0.0067	ug/L	1	-	U	Yes
gamma-Chlordane	0.0067	ug/L	1	-	U	Yes
Dieldrin	0.0067	ug/L	1	-	U	Yes
4,4'-DDD	0.0067	ug/L	1	-	U	Yes
4,4'-DDE	0.0067	ug/L	1	-	U	Yes
4,4'-DDT	0.0067	ug/L	1	-	U	Yes
Endrin	0.0067	ug/L	1	-	U	Yes
Endosulfan sulfate	0.0067	ug/L	1	-	U	Yes
Endrin aldehyde	0.0067	ug/L	1	-	U	Yes
Endrin ketone	0.0067	ug/L	1	-	U	Yes
Endosulfan-I	0.0067	ug/L	1	-	U	Yes
Endosulfan-II	0.0067	ug/L	1	-	U	Yes
Heptachlor	0.0067	ug/L	1	-	U	Yes
Heptachlor epoxide	0.0067	ug/L	1	-	U	Yes
Methoxychlor	0.013	ug/L	1	-	U	Yes
Toxaphene	0.17	ug/L	1	-	U	Yes

Sample ID: JC19023-4
Sample location: BMSMC Building 5 Area
Sampling date: 25-Apr-16
Matrix: Soil

METHOD: 8081B

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Aldrin	0.84	ug/kg	1	-	U	Yes
alpha-BHC	0.84	ug/kg	1	-	U	Yes
beta-BHC	0.84	ug/kg	1	-	U	Yes
delta-BHC	0.84	ug/kg	1	-	U	Yes
gamma-BHC (Lindane)	0.84	ug/kg	1	-	U	Yes
alpha-Chlordane	0.84	ug/kg	1	-	U	Yes
gamma-Chlordane	0.84	ug/kg	1	-	U	Yes
Dieldrin	0.84	ug/kg	1	-	U	Yes
4,4'-DDD	0.84	ug/kg	1	-	U	Yes
4,4'-DDE	0.84	ug/kg	1	-	U	Yes
4,4'-DDT	0.84	ug/kg	1	-	U	Yes
Endrin	0.84	ug/kg	1	-	U	Yes
Endosulfan sulfate	0.84	ug/kg	1	-	U	Yes
Endrin aldehyde	0.84	ug/kg	1	-	U	Yes
Endosulfan-I	0.84	ug/kg	1	-	U	Yes
Endosulfan-II	0.84	ug/kg	1	-	U	Yes
Heptachlor	0.84	ug/kg	1	-	U	Yes
Heptachlor epoxide	0.84	ug/kg	1	-	U	Yes
Methoxychlor	1.7	ug/kg	1	-	U	Yes
Endrin ketone	0.84	ug/kg	1	-	U	Yes
Toxaphene	21	ug/kg	1	-	U	Yes

DATA REVIEW WORKSHEETS

Project/Case Number: JC19023
 Sampling Date: April 22-25, 2016
 Shipping Date: April 25, 2016
 EPA Region No.: 2

REVIEW OF PESTICIDE ORGANIC PACKAGE

The following guidelines for evaluating volatile organics were created to delineate required validation actions. This document will assist the reviewer in using professional judgment to make more informed decision and in better serving the needs of the data users. The sample results were assessed according to USEPA data validation guidance documents in the following order of precedence *Hazardous Waste Support Section SOP No. HW-36A, Revision 0, June, 2015. SOM02.2. Pesticide Data Validation*. The QC criteria and data validation actions listed on the data review worksheets are from the primary guidance document, unless otherwise noted.

The hardcopied (laboratory name) Accutest data package received has been reviewed and the quality control and performance data summarized. The data review for VOCs included:

Lab. Project/SDG No.: JC19023 Sample matrix: Soil
 No. of Samples: 4
 Trip blank No.: -
 Field blank No.: -
 Equipment blank No.: JC19023-3
 Field duplicate No.: -
 Field spikes No.: -
 QC audit samples: -

<input checked="" type="checkbox"/> Data Completeness	<input checked="" type="checkbox"/> Laboratory Control Spikes
<input checked="" type="checkbox"/> Holding Times	<input checked="" type="checkbox"/> Field Duplicates
<input type="checkbox"/> N/A GC/MS Tuning	<input checked="" type="checkbox"/> Calibrations
<input checked="" type="checkbox"/> Internal Standard Performance	<input checked="" type="checkbox"/> Compound Identifications
<input checked="" type="checkbox"/> Blanks	<input checked="" type="checkbox"/> Compound Quantitation
<input checked="" type="checkbox"/> Surrogate Recoveries	<input checked="" type="checkbox"/> Quantitation Limits
<input checked="" type="checkbox"/> Matrix Spike/Matrix Spike Duplicate	

Overall Comments: TCL_pesticides_list_by_SW846-8081B

Definition of Qualifiers:

J- Estimated results
 U- Compound not detected
 R- Rejected data
 UJ- Estimated nondetect

Reviewer: Rafael Defaut
 Date: May 16, 2016

DATA REVIEW WORKSHEETS

DATA COMPLETENESS

MISSING INFORMATION

DATE LAB. CONTACTED

DATE RECEIVED

This image shows a blank sheet of white paper with horizontal ruling lines. A single diagonal line runs from the top-left corner towards the bottom-right corner, creating a margin. The lines are evenly spaced and extend across the width of the page. There is no handwriting or other markings on the paper.

DATA REVIEW WORKSHEETS

- e. Use professional judgment to qualify samples whose temperature upon receipt at the laboratory is either below 2 degrees centigrade or above 6 degrees centigrade.
- f. If technical holding times are grossly exceeded, use professional judgment to qualify the data.

Qualify non-aqueous sample results using preservation and technical holding time information as follows:

- a. If there is no evidence that the samples were properly preserved ($T = 4^{\circ}\text{C} \pm 2^{\circ}\text{C}$), and the samples were extracted or analyzed within the technical holding time, qualify detects as estimated (J) and non-detects as estimated (UJ).
- b. If there is no evidence that the samples were properly preserved ($T = 4^{\circ}\text{C} \pm 2^{\circ}\text{C}$), and the samples were extracted or analyzed outside the technical holding time, qualify detects as estimated (J) and non-detects as estimated (UJ).
- c. If the samples were properly preserved, and were extracted and analyzed within the technical holding time, no qualification of the data is necessary.
- d. If the samples were properly preserved, and were extracted or analyzed outside the technical holding time, qualify detects as estimated (J) and non-detects as estimated (UJ). Note in the Data Review Narrative that holding times were exceeded and the effect of exceeding the holding time on the resulting data.
- e. Use professional judgment to qualify samples whose temperature upon receipt at the laboratory is either below 2 degrees centigrade or above 6 degrees centigrade.
- f. If technical holding times are grossly exceeded, use professional judgment to qualify the data.

DATA REVIEW WORKSHEETS

All criteria were met X
Criteria were not met see below

GAS CHROMATOGRAPH WITH ELECTRON CAPTURE DETECTOR (GC/ECD) INSTRUMENT PERFORMANCE CHECK (SECTIONS 1 TO 5)

1. Resolution Check Mixture

Criteria

Is the resolution between two adjacent peaks in the Resolution Check Mixture C greater than or equal to 80.0% for all analytes for the primary column and greater than or equal to 50.0% for the confirmation column? Yes? or No?

Is the resolution between two adjacent peaks in the Resolution Check Mixture (A and B) greater than or equal to 60.0%? Yes? or No?

Note: If resolution criteria are not met, the quantitative results may not be accurate due to inadequate resolution. Qualitative identifications may also be questionable if coelution exists.

Action

- a. Qualify detects for target compounds that were not adequately resolved as tentatively identified (NJ).
- b. Qualify non-detected compounds as unusable (R).

2. Performance Evaluation Mixture (PEM) Resolution Criteria

Criteria

Is PEM analysis performed at the required frequency (at the end of each pesticide initial calibration sequence and every 12 hours)? Yes? or No?

Action

- a. If PEM is not performed at the required frequency, qualify all associated sample and blank results as unusable (R).

Criteria

Is PEM % Resolution < 90%? Yes? or No?

Action

- a. a. Qualify detects for target compounds that were not adequately resolved as tentatively identified (NJ).
- b. Qualify non-detected compounds as unusable (R).

All criteria were met ____ N/A ____
Criteria were not met see below ____

3. PEM 4,4'-DDT Breakdown

Criteria

Is the PEM 4,4'-DDT % Breakdown >20.0% and 4,4'-DDT is detected? Yes? or No?

Action

a. Qualify detects for 4,4'-DDT; detects for 4,4'-DDD; and detects for 4,4'-DDE as estimated (J)

Criteria

Is the PEM 4,4'-DDT % Breakdown >20.0% and 4,4'-DDT is not detected Yes? or No?

Action

- a. Qualify non-detects for 4,4'- DDT as unusable (R)
- b. Qualify detects for 4,4'-DDD as tentatively identified (NJ)
- c. Qualify detects for 4,4'-DDE as tentatively identified (NJ)

4. PEM Endrin Breakdown

Criteria

Is the PEM Endrin % Breakdown >20.0% and Endrin is detected? Yes? or No?

Action

a. Qualify detects for Endrin; detects for Endrin aldehyde; and detects for Endrin ketone as estimated (J)

Criteria

Is the PEM Endrin % Breakdown >20.0% and Endrin is not detected Yes? or No?

Action

- a. Qualify non-detects for Endrin as unusable (R)
- b. Qualify detects for Endrin aldehyde as tentatively identified (NJ)
- c. Qualify detects for Endrin ketone as tentatively identified (NJ)

All criteria were met X
Criteria were not met see below

5. Mid-point Individual Standard Mixture Resolution -

Criteria

Is the resolution between two adjacent peaks in the Resolution Check Mixture C greater than or equal to 80.0% for all analytes for the primary column and greater than or equal to 50.0% for the confirmation column?
Yes? or No?

Is the resolution between two adjacent peaks in the Resolution Check Mixture (A and B) greater than or equal to 90.0%?
Yes? or No?

Note: If resolution criteria are not met, the quantitative results may not be accurate due to inadequate resolution. Qualitative identifications may also be questionable if coelution exists.

Action

- a. Qualify detects for target compounds that were not adequately resolved as tentatively identified (NJ).
- b. Qualify non-detected compounds as unusable (R).

Criteria

Is mid-point individual standard mixture analysis performed at the required frequency (every 12 hours)?
Yes? or No?

Action

- a. If the mid-point individual standard mixture analysis is not performed at the required frequency, qualify all associated sample and blank results as unusable (R).

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

CALIBRATION VERIFICATION

Compliance requirements for satisfactory instrument calibration are established to ensure that the instrument is capable of producing and maintaining acceptable quantitative data.

Date of initial calibration: 04/27/16 04/18/16
 Dates of initial calibration verification: 04/27/16 04/18/16
 Dates of continuing calibration: 04/29/16 04/27/16
 Dates of final calibration: 04/29/16 04/27/16
 Instrument ID numbers: GC1G HP_G153A
 Matrix/Level: Aqueous/low Aqueous/low

DATE	LAB ID#	FILE	CRITERIA OUT RFs, %RSD, %D, r	COMPOUND	SAMPLES AFFECTED
Initial and initial calibration verification within the guidance document performance criteria. Continuing calibration % differences meet the performance criteria.					

Criteria

Are a five point calibration curve delivered with concentration levels as shown in Table 3 of SOP HW-36A, Revision 0, June, 2015? Yes? or No?

Actions

If the standard concentrations listed in Table 3 are not used, use professional judgment to evaluate the effect on the data

Criteria

Are RT Windows calculated correctly? Yes? or No?

Action

Recalculate the windows and use the corrected values for all evaluations.

Criteria

Are the Percent Relative Standard Deviation (%RSD) of the CFs for each of the single component target compounds less than or equal to 20.0%, except for alpha-BHC and delta-BHC?

Yes? or No?

Are the %RSD of the CFs for alpha-BHC and delta-BHC less than or equal to 25.0%. Yes? or No?

DATA REVIEW WORKSHEETS

Is the %RSD of the CFs for each of the Toxaphene peaks must be < 30% when 5-point ICAL is performed? **Yes?** or **No?**

Is the %RSD of the CFs for the two surrogates (tetrachloro-m-xylene and decachlorobiphenyl) less than or equal to 30.0%. **Yes?** or **No?**

Action

- a. If the %RSD criteria are not met, qualify detects as estimated (J) and use professional judgment to qualify non-detected target compounds.
- b. If the %RSD criteria are within allowable limits, no qualification of the data is necessary

Continuing Calibration Checks

Criteria

Is the continuing calibration standard analyzed at the acceptable time intervals? **Yes?** or **No?**

Action

- a. If more than 14 hours has elapsed from the injection of the instrument blank that begins an analytical sequence (opening CCV) and the injection of either a PEM or mid-point concentration of the Individual Standard Mixtures (A and B) or (C), qualify all data as unusable (R).
- b. If more than 12 hours has elapsed from the injection of the instrument blank that begins an analytical sequence (opening CCV) and the injection of the last sample or blank that is part of the same analytical sequence, qualify all data as unusable (R).
- c. If more than 72 hours has elapsed from the injection of the sample with a Toxaphene detection and the Toxaphene Calibration Verification Standard (CS3), qualify all data as unusable (R).

Criteria

Is the Percent Difference (%D) within $\pm 25.0\%$ for the PEM sample? **Yes?** or **No?**

Action

- a. Qualify associated detects as estimated (J) and non-detects as estimated (UJ).

Criteria

For the Calibration Verification Standard (CS3); is the Percent Difference (%D) within $\pm 25.0\%$? **Yes?** or **No?**

Action

Qualify associated detects as estimated (J) and non-detects as estimated (UJ).

DATA REVIEW WORKSHEETS

All criteria were met N/A
Criteria were not met
and/or see below

Criteria

Is the PEM 4,4'-DDT % Breakdown >20.0% and 4,4'-DDT is detected? Yes? or No?

Action

- a. Qualify detects for 4,4'-DDT; detects for 4,4'-DDD; and detects for 4,4'-DDE as estimated (J)
- b. Non-detected associated compounds are not qualified

Criteria

Is the PEM 4,4'-DDT % Breakdown >20.0% and 4,4'-DDT is not detected Yes? or No?

Action

- a. Qualify non-detects for 4,4'- DDT as unusable (R)
- b. Qualify detects for 4,4'-DDD as tentatively identified (NJ)
- c. Qualify detects for 4,4'-DDE as tentatively identified (NJ)

Criteria

Is the PEM Endrin % Breakdown >20.0% and Endrin is detected? Yes? or No?

Action

- a. Qualify detects for Endrin; detects for Endrin aldehyde; and detects for Endrin ketone as estimated (J)
- b. Non-detected associated compounds are not qualified

Criteria

Is the PEM Endrin % Breakdown >20.0% and Endrin is not detected Yes? or No?

Action

- a. Qualify non-detects for Endrin as unusable (R)
- b. Qualify detects for Endrin aldehyde as tentatively identified (NJ)
- c. Qualify detects for Endrin ketone as tentatively identified (NJ)

A separate worksheet should be filled for each initial curve

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

BLANK ANALYSIS RESULTS (Sections 1 & 2)

The assessment of the blank analysis results is to determine the existence and magnitude of contamination problems. The criteria for evaluation of blanks apply only to blanks associated with the samples, including trip, equipment, and laboratory blanks. If problems with any blanks exist, all data associated with the case must be carefully evaluated to determine whether or not there is an inherent variability in the data for the case, or if the problem is an isolated occurrence not affecting other data.

List the contamination in the blanks below. High and low levels blanks must be treated separately.

CRQL concentration N/A

Laboratory blanks

DATE ANALYZED	LAB ID	LEVEL/ MATRIX	COMPOUND	CONCENTRATION UNITS
<u>No target analytes detected in method blanks at a reporting limit of 0.01 and 0.001 ug/L.</u>				

Field/Equipment/Trip blank

DATE ANALYZED	LAB ID	LEVEL/ MATRIX	COMPOUND	CONCENTRATION UNITS
<u>No field/trip blanks analyzed with this data package. No target analytes detected in the equipment blank.</u>				

DATA REVIEW WORKSHEETS

All criteria were met ☒ X
 Criteria were not met
 and/or see below _____

BLANK ANALYSIS RESULTS (Section 3)

Blank Actions

Action Levels (ALs) should be based upon the highest concentration of contaminant determined in any blank. Do not qualify any blank with another blank. The ALs for samples which have been diluted should be corrected for the sample dilution factor and/or % moisture, where applicable. No positive sample results should be reported unless the concentration of the compound in the samples exceeds the ALs:

The concentration of non-target compounds in all blanks must be less than or equal to 10 µg/L. The concentration of each target compound found in the method or field blanks must be less than its CRQL listed in the method.

Data concerning the field blanks are not evaluated as part of the CCS process. If field blanks are present, the data reviewer should evaluate this data in a similar fashion as the method blanks.

Specific actions are as follows:

Blank Actions for Pesticide Analyses

Blank Type	Blank Result	Sample Result	Action for Samples
Method, Sulfur Cleanup, Instrument, Field, TCLP/SPLP	Detects	Not detected	No qualification required
	< CRQL	< CRQL	Report CRQL value with a U
		≥ CRQL	No qualification required
	> CRQL	< CRQL	Report CRQL value with a U
		≥ CRQL and ≤ blank concentration	Report blank value for sample concentration with a U
		≥ CRQL and > blank concentration	No qualification required
	= CRQL	≤ CRQL	Report CRQL value with a U
		> CRQL	No qualification required
	Gross contamination	Detects	Report blank value for sample concentration with a U

DATA REVIEW WORKSHEETS

All criteria were met ☒ X
 Criteria were not met
 and/or see below _____

SURROGATE SPIKE RECOVERIES

Laboratory performance of individual samples is established by evaluation of surrogate spike recoveries. All samples are spiked with surrogate compounds prior to sample analysis. The accuracy of the analysis is measured by the surrogate percent recovery. Since the effects of the sample matrix are frequently outside the control of the laboratory and may present relatively unique problems, the validation of data is frequently subjective and demands analytical experience and professional judgment.

List the percent recoveries (%Rs) which do not meet the criteria for surrogate recovery.

Matrix: Aqueous

Lab Sample ID	Lab File ID	S1 a	S1 b	S2 a	S2 b
JC19023-3	6G34452.D	96	97	61	67
OP93399-BS1	6G34448.D	66	68	68	73
OP93399-BSD	6G34449.D	93	95	94	100
OP93399-MB1	6G34447.D	92	99	98	103

Surrogate Compounds Recovery Limits

S1 = Tetrachloro-m-xylene 26-132%
 S2 = Decachlorobiphenyl 10-118%

(a) Recovery from GC signal #1

(b) Recovery from GC signal #2

Matrix: Soil

Lab Sample ID	Lab File ID	S1 a	S1 b	S2 a	S2 b
JC19023-1	1G122627.D	59	65	59	69
JC19023-2	1G122628.D	80	84	76	82
JC19023-4	1G122629.D	61	58	58	57
OP93471-BS1	1G122605.D	71	74	75	76
OP93471-MB1	1G122604.D	66	69	70	71
OP93471-MS	1G122619.D	60	62	52	57
OP93471-MSD	1G122620.D	80	82	68	76

Surrogate Compounds Recovery Limits

S1 = Tetrachloro-m-xylene 24-136%
 S2 = Decachlorobiphenyl 10-153%

(a) Recovery from GC signal #1

(b) Recovery from GC signal #2

Note: Surrogate recoveries within laboratory control limits.

DATA REVIEW WORKSHEETS

Actions:

- a. For any surrogate recovery greater than 150%, qualify detected target compounds as biased high (J+).
- b. Do not qualify non-detected target compounds for surrogate recovery > 150 %.
- c. If both surrogate recoveries are greater than or equal to 30% and less than or equal to 150%, no qualification of the data is necessary.
- d. For any surrogate recovery greater than or equal to 10% and less than 30%, qualify detected target compounds as biased low (J-).
- e. For any surrogate recovery greater than or equal to 10% and less than 30%, qualify non-detected target compounds as approximated (UJ).
- f. If low surrogate recoveries are from sample dilution, professional judgment should be used to determine if the resulting data should be qualified. If sample dilution is not a factor:
 - i. Qualify detected target compounds as biased low (J-).
 - ii. Qualify non-detected target compounds as unusable (R).
- g. If surrogate RTs in PEMs, Individual Standard Mixtures, samples, and blanks are outside of the RT Windows, the reviewer must use professional judgment to qualify data.
- h. If surrogate RTs are within RT windows, no qualification of the data is necessary.
- i. If the two surrogates were not added to all samples, MS/MSDs, standards, LCSs, and blanks, use professional judgment in qualifying data as missing surrogate analyte may not directly apply to target analytes.

Summary Surrogate Actions for Pesticide Analyses

Criteria	Action*	
	Detected Target Compounds	Non-detected Target Compounds
%R > 150%	J+	No qualification
30% < %R < 150%	No qualification	
10% < %R < 30%	J-	UJ
%R < 10% (sample dilution not a factor)	J-	R
%R < 10% (sample dilution is a factor)	Use professional judgment	
RT out of RT window	Use professional judgment	
RT within RT window	No qualification	

- * Use professional judgment in qualifying data, as surrogate recovery problems may not directly apply to target analytes.

All criteria were met ___N/A___
Criteria were not met
and/or see below _____

MATRIX SPIKE/MATRIX SPIKE DUPLICATE (MS/MSD)

This data is generated to determine long term precision and accuracy in the analytical method for various matrices. This data alone cannot be used to evaluate the precision and accuracy of individual samples. If any % R in the MS or MSD falls outside the designated range, the reviewer should determine if there are matrix effects, i.e. LCS data are within the QC limits but MS/MSD data are outside QC limit.

1. MS/MSD Recoveries and Precision Criteria

Data for MS and MSDs will not be present unless requested by the Region.

Notify the Contract Laboratory Program Project Officer (CLP PO) if a field blank was used for the MS and MSD, unless designated as such by the Region.

NOTE: For a Matrix Spike that does not meet criteria, apply the action to only the field sample used to prepare the Matrix Spike sample. If it is clearly stated in the data validation materials that the samples were taken through incremental sampling or some other method guaranteeing the homogeneity of the sample group, then the entire sample group may be qualified.

List the %Rs, RPD of the compounds which do not meet the criteria.

Sample ID: ___JC19164-5_____

Matrix/Level: ___Soil_____

MS OR MSD	COMPOUND	% R	RPD	QC LIMITS	ACTION
-----------	----------	-----	-----	-----------	--------

___MS/MSD_%_recoveries_and_RPD_within_laboratory_control_limits._____

Note: No MS/MSD analyzed with this data package for aqueous matrix. Blank spike/blank spike duplicate used to assess accuracy. % recoveries and RPD within laboratory control limits.

Action

No qualification of the data is necessary on MS and MSD data alone. However, using professional judgment, the validator may use the MS and MSD results in conjunction with other QC criteria and determine the need for some qualification of the data.

A separate worksheet should be used for each MS/MSD pair.

All criteria were met X
 Criteria were not met
 and/or see below

LABORATORY CONTROL SAMPLE (LCS) ANALYSIS

This data is generated to determine accuracy of the analytical method for various matrices.

1. LCS Recoveries Criteria

LCS Spike Compound	Recovery Limits (%)
gamma-BHC	50 – 120
Heptachlor epoxide	50 – 150
Dieldrin	30 – 130
4,4'-DDE	50 – 150
Endrin	50 – 120
Endosulfan sulfate	50 – 120
trans-Chlordane	30 – 130
Tetrachloro-m-xylene (surrogate)	30 – 150
Decachlorobiphenyl (surrogate)	30 – 150

LCS concentrations: 0.167 ug/L

List the %R of compounds which do not meet the criteria

LCS ID	COMPOUND	% R	QC LIMIT

Action

The following guidance is suggested for qualifying sample data for which the associated LCS does not meet the required criteria.

- If the LCS recovery exceeds the upper acceptance limit, qualify detected target compounds as estimated (J). Do not qualify non-detected target compounds.
- If the LCS recovery is less than the lower acceptance limit, qualify detected target compounds as estimated (J) and non-detects as unusable (R).
- Use professional judgment to qualify data for compounds other than those compounds that are included in the LCS.
- Use professional judgment to qualify non-LCS compounds. Take into account the compound class, compound recovery efficiency, analytical problems associated with each compound, and comparability in the performance of the LCS compound to the non-LCS compound.
- If the LCS recovery is within allowable limits, no qualification of the data is necessary.

DATA REVIEW WORKSHEETS

2. Frequency Criteria:

Where LCS analyzed at the required frequency and for each matrix? **Yes** or No.

If no, the data may be affected. Use professional judgment to determine the severity of the effect and qualify data accordingly. Discuss any actions below and list the samples affected.

Note: Analytes do not meet the RPD in the blank spike/blank spike duplicate (aqueous matrix), see below. No action taken, professional judgment.

The QC reported here applies to the following samples:
JC19023-3

Method: **SW846 8081B**

Compound	Spike ug/l	BSP ug/l	BSP %	BSD ug/l	BSD %	RPD	Limits Rec/RPD
Aldrin	0.167	0.099	59	0.15	90	41* a	38-138/30
alpha-BHC	0.167	0.10	60	0.15	90	40* a	51-144/30
beta-BHC	0.167	0.11	66	0.15	90	31* a	55-137/30
delta-BHC	0.167	0.11	66	0.16	96	37* a	33-148/30
gamma-BHC (Lindane)	0.167	0.11	66	0.15	90	31* a	53-142/30
alpha-Chlordane	0.167	0.10	60	0.15	90	40* a	52-137/30
gamma-Chlordane	0.167	0.098	59	0.15	90	42* a	52-140/30
Dieldrin	0.167	0.10	60	0.15	90	40* a	54-144/30
4,4'-DDD	0.167	0.11	66	0.15	90	31* a	46-151/30
4,4'-DDE	0.167	0.10	60	0.15	90	40* a	51-142/30
4,4'-DDT	0.167	0.11	66	0.16	96	37* a	53-158/30
Endrin	0.167	0.11	66	0.16	96	37* a	53-149/30
Endosulfan sulfate	0.167	0.12	72	0.17	102	34* a	49-148/30
Endrin aldehyde	0.167	0.11	66	0.15	90	31* a	45-137/30
Endrin ketone	0.167	0.11	66	0.16	96	37* a	48-142/30
Endosulfan-I	0.167	0.10	60	0.15	90	40* a	54-139/30
Endosulfan-II	0.167	0.11	66	0.15	90	31* a	54-142/30
Heptachlor	0.167	0.10	60	0.15	90	40* a	45-137/30
Heptachlor epoxide	0.167	0.11	66	0.15	90	31* a	54-140/30
Methoxychlor	0.167	0.12	72	0.17	102	34* a	51-157/30

(a) Analytical precision exceeds in-house control limits.

DATA REVIEW WORKSHEETS

All criteria were met _____
Criteria were not met _____
and/or see below N/A

FLORISIL CARTRIDGE PERFORMANCE CHECK

NOTE: Florisil cartridge cleanup is mandatory for all extracts.

Criteria

Is the Florisil cartridge performance check conducted at least once on each lot of cartridges used for sample cleanup or every 6 months, whichever is most frequent? Yes? or No? N/A

Criteria

Are the results for the Florisil Cartridge Performance Check solution included with the data package? Yes? or No? N/A

Note: If % criteria are not met, examine the raw data for the presence of polar interferences and use professional judgment in qualifying the data as follows:

Action:

- a. If the Percent Recovery is greater than 120% for any of the pesticide target compounds in the Florisil Cartridge Performance Check, qualify detected compounds as estimated (J). Do not qualify non-detected target compounds.
- b. If the Percent Recovery is greater than or equal to 80% and less than or equal to 120% for all the pesticide target compounds, no qualification of the data is necessary.
- c. If the Percent Recovery is greater than or equal to 10% and less than 80% for any of the pesticide target compounds in the Florisil Cartridge Performance Check, qualify detected target compounds as estimated (J) and non-detected target compounds as approximated (UJ).
- d. If the Percent Recovery is less than 10% for any of the pesticide target compounds in the Florisil Cartridge Performance Check, qualify detected compounds as estimated (J) and qualify non-detected target compounds as unusable (R).
- e. If the Percent Recovery of 2,4,5-trichlorophenol in the Florisil Cartridge Performance Check is greater than or equal to 5%, use professional judgment to qualify detected and non-detected target compounds, considering interference on the sample chromatogram.

Note: State in the Data Review Narrative potential effects on the sample data resulting from the Florisil Cartridge Performance Check analysis not yielding acceptable results.

Note: No information for florisil cartridge performance check included in data package. Florisil cartridge used for sample extraction/clean-up. No qualification of the data performed, professional judgment.

All criteria were met N/A
Criteria were not met
and/or see below

GEL PERMEATION CHROMATOGRAPHY (GPC) PERFORMANCE CHECK

NOTE: GPC cleanup is mandatory for all soil samples.

If GPC criteria are not met, examine the raw data for the presence of high molecular weight contaminants; examine subsequent sample data for unusual peaks; and use professional judgment in qualifying the data. Notify the Contract Laboratory Program Project Officer (CLP PO) if the laboratory chooses to analyze samples under unacceptable GPC criteria.

Action:

- a. If the Percent Recovery is less than 10% for the pesticide compounds and surrogates during the GPC calibration check, the non-detected target compounds may be suspect, qualify detected compounds as estimated (J).
- b. If the Percent Recovery is less than 10% for the pesticide compounds and surrogates during the GPC calibration check, qualify all non-detected target compounds as unusable (R).
- c. If the Percent Recovery is greater than or equal to 10% and is less than 80% for any of the pesticide target compounds in the GPC calibration, qualify detected target compounds as estimated (J) and non-detected target compounds as approximated (UJ).
- d. If the Percent Recovery is greater than or equal to 80% and less than or equal to 120% for all the pesticide target compounds, no qualification of the data is necessary.
- e. If high recoveries (i.e., greater than 120%) were obtained for the pesticides and surrogates during the GPC calibration check, qualify detected compounds as estimated (J). Do not qualify non-detected target compounds.

Note: State in the Data Review Narrative potential effects on the sample data resulting from the GPC cleanup analyses not yielding acceptable results.

Note: No information for performance of GPC cleanup included in data package. No qualification of the data performed, professional judgment.

DATA REVIEW WORKSHEETS

All criteria were met X
Criteria were not met
and/or see below _____

TARGET COMPOUND IDENTIFICATION

Criteria:

1. Is Retention Times (RTs) of both of the surrogates and reported target compounds in each sample within the calculated RT Windows on both columns? Yes? or No?
2. Is the Tetrachloro-m-xylene (TCX) RT ± 0.05 minutes of the Mean RT (RT) determined from the initial calibration and Decachlorobiphenyl (DCB) within ± 0.10 minutes of the RT determined from the initial calibration? Yes? or No?
3. Is the Percent Difference (%D) for the detected mean concentrations of a pesticide target compound between the two Gas Chromatograph (GC) columns within the inclusive range of ± 25.0 %? Yes? or No?
4. When no analytes are identified in a sample; are the chromatograms from the analyses of the sample extract and the low-point standard of the initial calibration associated with those analyses on the same scaling factor? Yes? or No?
5. Does the chromatograms display the Single Component Pesticides (SCPs) detected in the sample and the largest peak of any multi-component analyte detected in the sample at less than full scale. Yes? or No?
6. If an extract is diluted; does the chromatogram display SCPs peaks between 10-100% of full scale, and multi-component analytes between 25-100% of full scale? Yes? or No? N/A
7. For any sample; does the baseline of the chromatogram return to below 50% of full scale before the elution time of alpha-BHC, and also return to below 25% of full scale after the elution time of alpha-BHC and before the elution time of DCB? Yes? or No?
8. If a chromatogram is replotted electronically to meet these requirements; is the scaling factor used displayed on the chromatogram, and both the initial chromatogram and the replotted chromatogram submitted in the data package. Yes? or No?

Action:

- a. If the qualitative criteria for both columns were not met, all target compounds that are reported as detected should be considered non-detected.
- b. Use professional judgment to assign an appropriate quantitation limit using the following guidance:
 - i. If the detected target compound peak was sufficiently outside the pesticide RT Window, the reported values may be a false positive and should be replaced with the sample Contract Required Quantitation Limits (CRQL) value.

DATA REVIEW WORKSHEETS

- ii. If the detected target compound peak poses an interference with potential detection of another target peak, the reported value should be considered and qualified as unusable (R).

c. If the data reviewer identifies a peak in both GC column analyses that falls within the appropriate RT Windows, but was reported as a non-detect, the compound may be a false negative. Use professional judgment to decide if the compound should be included.

Note: State in the Data Review Narrative all conclusions made regarding target compound identification.

d. If the Toxaphene peak RT windows determined from the calibration overlap with SCPs or chromatographic interferences, use professional judgment to qualify the data.

e. If target compounds were detected on both GC columns, and the Percent Difference between the two results is greater than 25.0%, consider the potential for coelution and use professional judgment to decide whether a much larger concentration obtained on one column versus the other indicates the presence of an interfering compound. If an interfering compound is indicated, use professional judgment to determine how best to report, and if necessary, qualify the data according to these guidelines.

f. If Toxaphene exhibits a marginal pattern-matching quality, use professional judgment to establish whether the differences are due to environmental "weathering" (i.e., degradation of the earlier eluting peaks relative to the later eluting peaks). If the presence of Toxaphene is strongly suggested, report results as presumptively present (N).

GAS CHROMATOGRAPH/MASS SPECTROMETER (GC/MS) CONFIRMATION

NOTE: This confirmation is not usually provided by the laboratory. In cases where it is provided, use professional judgment to determine if data qualified with "C" can be salvaged if it was previously qualified as unusable (R).

Action:

a. If the quantitative criteria for both columns were met (≥ 5.0 ng/ μ L for SCPs and ≥ 125 ng/ μ L for Toxaphene), determine whether GC/MS confirmation was performed. If it was performed, qualify the data using the following guidance:

- i. If GC/MS confirmation was not required because the quantitative criteria for both columns was not met, but it was still performed, use professional judgment when evaluating the data to decide whether the detect should be qualified with "C".
- ii. If GC/MS confirmation was performed, but unsuccessful for a target compound detected by GC/ECD analysis, qualify those detects as "X".

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

COMPOUND QUANTITATION AND REPORTED CONTRACT REQUIRED QUANTITATION LIMITS (CRQLS)

The sample quantitation evaluation is to verify laboratory quantitation results. In the space below, please show a minimum of one sample calculation:

op93471-ms (Matrix Spike) alpha-BHC RF = 1.286

$$\begin{aligned} [] &= (23933340)(50)/(66473745)(1.286) \\ &= 13.97 \text{ ppb} \quad \text{Ok} \end{aligned}$$

Action:

- If sample quantitation is different from the reported value, qualify result as unusable (R).
- When a sample is analyzed at more than one dilution, the lowest CRQLs are used unless a QC exceedance dictates the use of the higher CRQLs from the diluted sample.
- Replace concentrations that exceed the calibration range in the original analysis by crossing out the "E" and its corresponding value on the original reporting form and substituting the data from the diluted sample.
- Results between the MDL and CRQL should be qualified as estimated (J).
- Results less than the MDL should be reported at the CRQL and qualified (U). MDLs themselves are not reported.
- For non-aqueous samples, if the percent moisture is less than 70.0%, no qualification of the data is necessary. If the percent moisture is greater than or equal to 70.0% and less than 90.0%, qualify detects as estimated (J) and non-detects as approximated (UJ). If the percent moisture is greater than or equal to 90.0%, qualify detects as estimated (J) and non-detects as unusable (R) (see Table).

Percent Moisture Actions for Pesticide Analysis for Non-Aqueous Samples

Criteria	Action	
	Detected Associated Compounds	Non-detected Associated Compounds
% Moisture < 70.0	No qualification	
70.0 < % Moisture < 90.0	J	UJ
% Moisture > 90.0	J	R

DATA REVIEW WORKSHEETS

List samples which have $\leq 50\%$ solids

Note: If any discrepancies are found, the Region's designated representative may contact the laboratory to obtain additional information that could resolve any differences. If a discrepancy remains unresolved, the reviewer must use professional judgment to decide which value is the most accurate. Under these circumstances, the reviewer may determine that qualification of data is warranted. Note in the Data Review Narrative a description of the reasons for data qualification and the qualification that is applied to the data.

Dilution performed

SAMPLE ID	DILUTION FACTOR	REASON FOR DILUTION

DATA REVIEW WORKSHEETS

All criteria were met N/A
 Criteria were not met
 and/or see below

FIELD DUPLICATE PRECISION

NOTE: In the absence of QAPP guidance for validating data from field duplicates, the following action will be taken.

Field duplicate samples may be taken and analyzed as an indication of overall precision. These analyses measure both field and lab precision; therefore, the results may have more variability than laboratory duplicates which only laboratory performance. It is also expected that soil duplicate results will have a greater variance than water matrices due to difficulties associated with collecting identical field duplicate samples. Identify which samples within the data package are field duplicates. Estimate the relative percent difference (RPD) between the values for each compound. If large RPDs (> 50%) is observed, confirm identification of samples and note difference in the executive summary.

Sample IDs: -

Matrix: -

COMPOUND	SQL ug/L	SAMPLE CONC.	DUPLICATE CONC.	RPD	ACTION
No field/laboratory duplicate analyzed with this data package. LCS/LCSD % recoveries RPD used to assess precision. RPD within the required criteria of < 50 %.					

Actions:

a. Qualify as estimated positive results (J) and nondetects (UJ) for the compound that exceeded the above criteria. For organics, only the sample and duplicate will be qualified.

b. If an RPD cannot be calculated because one or both of the sample results is not detected, the following actions apply:

- i. If one sample result is not detected and the other is greater than 5x the SQL qualify (J/UJ).
- ii. If one sample value is not detected and the other is greater than 5x the SQL and the SQLs for the sample and duplicate are significantly different, use professional judgment to determine if qualification is appropriate.
- iii. If one sample value is not detected and the other is less than 5x, use professional judgment to determine if qualification is appropriate.
- iv. If both sample and duplicate results are not detected, no action is needed.

DATA REVIEW WORKSHEETS

OVERALL ASSESSMENT OF DATA

Action:

1. Use professional judgment to determine if there is any need to qualify data which were not qualified based on the Quality Control (QC) criteria previously discussed.
2. Write a brief narrative to give the user an indication of the analytical limitations of the data.

Note: The Contract Laboratory Program Project Officer (CLP PO) must be informed if any inconsistency of the data with the Sample Delivery Group (SDG) Narrative. If sufficient information on the intended use and required quality of the data is available, the reviewer should include their assessment of the usability of the data within the given context. This may be used as part of a formal Data Quality Assessment (DQA).

Overall assessment of the data: Results are valid; the data can be used for decision making purposes.

EXECUTIVE NARRATIVE

SDG No: **JC19023** Laboratory: **Accutest, Florida**
Analysis: **SW846-8015C** Number of Samples: **4**
Location: **BMSMC, Building 5 Area**
Humacao, PR

SUMMARY: Three (3) soil samples and one (1) equipment blank sample were analyzed for the low molecular weight alcohols (LMWAs) list following method SW846-8015C. The sample results were assessed according to USEPA data validation guidance documents in the following order of precedence: "Test Methods for Evaluating Solid Waste, Physical/Chemical Methods SW-846 (Final Update III, December 1996)," specifically for Methods 8000/8015C are utilized. The QC criteria and data validation actions listed on the data review worksheets are from the primary guidance document, unless otherwise noted.

Results are valid and can be used for decision making purposes.


Critical issues: **None**
Major: **None**
Minor: **None**

Critical findings: **None**
Major findings: **None**
Minor findings: 1. All samples analyzed within the recommended method holding time. All samples properly preserved except samples JC19023-1; JC19023-2; and JC18649-4 (Soil) that were not preserved within 48 hours of sampling. Results qualified as estimated (UJ) in the affected sample.
2. JC18972-2A MS/MSD recoveries outside the laboratory control limits but within generally acceptable control limits. No action taken, qualification applies only to the unspiked sample. Sample JC18972-2A from another data package, used for QC purposes only.

COMMENTS: Results are valid and can be used for decision making purposes.

Reviewers Name: **Rafael Infante**
Chemist License 1888

Signature:

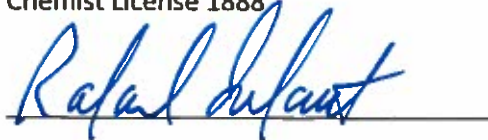

Date: **May 16, 2016**

SAMPLE ORGANIC DATA SAMPLE SUMMARY**Sample ID: JC19023-1A****Sample location: BMSMC Building 5 Area****Sampling date: 4/22/2016****Matrix: Soil****METHOD: 8015C**

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Ethanol	13	mg/kg	1.0	-	UJ	Yes
Isobutyl Alcohol	13	mg/kg	1.0	-	UJ	Yes
Isopropyl Alcohol	13	mg/kg	1.0	-	UJ	Yes
n-Propyl Alcohol	13	mg/kg	1.0	-	UJ	Yes
n-Butyl Alcohol	13	mg/kg	1.0	-	UJ	Yes
Methanol	13	mg/kg	1.0	-	UJ	Yes

Sample ID: JC19023-2A**Sample location: BMSMC Building 5 Area****Sampling date: 4/22/2016****Matrix: Soil****METHOD: 8015C**

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Ethanol	13	mg/kg	1.0	-	UJ	Yes
Isobutyl Alcohol	13	mg/kg	1.0	-	UJ	Yes
Isopropyl Alcohol	13	mg/kg	1.0	-	UJ	Yes
n-Propyl Alcohol	13	mg/kg	1.0	-	UJ	Yes
n-Butyl Alcohol	13	mg/kg	1.0	-	UJ	Yes
Methanol	13	mg/kg	1.0	-	UJ	Yes

Sample ID: JC19023-3**Sample location: BMSMC Building 5 Area****Sampling date: 4/22/2016****Matrix: AQ - Equipment Blank****METHOD: 8015C**

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Ethanol	5.0	mg/l	1.0	-	U	Yes
Isobutyl Alcohol	5.0	mg/l	1.0	-	U	Yes
Isopropyl Alcohol	5.0	mg/l	1.0	-	U	Yes
n-Propyl Alcohol	5.0	mg/l	1.0	-	U	Yes
n-Butyl Alcohol	5.0	mg/l	1.0	-	U	Yes
Methanol	5.0	mg/l	1.0	-	U	Yes

Sample ID: JC19023-4
Sample location: BMSMC Building 5 Area
Sampling date: 4/5/2016
Matrix: Soil

METHOD: 8015C

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Ethanol	13	mg/kg	1.0	-	UJ	Yes
Isobutyl Alcohol	13	mg/kg	1.0	-	UJ	Yes
Isopropyl Alcohol	13	mg/kg	1.0	-	UJ	Yes
n-Propyl Alcohol	13	mg/kg	1.0	-	UJ	Yes
n-Butyl Alcohol	13	mg/kg	1.0	-	UJ	Yes
Methanol	13	mg/kg	1.0	-	UJ	Yes

DATA REVIEW WORKSHEETS

Project Number: JC19023

Date: 04/22-25/2016

Shipping Date: 04/25/2016

EPA Region: 2

REVIEW OF VOLATILE ORGANIC PACKAGE

The following guidelines for evaluating volatile organics were created to delineate required validation actions. This document will assist the reviewer in using professional judgment to make more informed decision and in better serving the needs of the data users. The sample results were assessed according to USEPA data validation guidance documents in the following order of precedence: "Test Methods for Evaluating Solid Waste, Physical/Chemical Methods SW-846 (Final Update III, December 1996)," specifically for Methods 8000/8015C are utilized. The QC criteria and data validation actions listed on the data review worksheets are from the primary guidance document, unless otherwise noted.

The hardcopied (laboratory name) Accutest data package received has been reviewed and the quality control and performance data summarized. The modified data review for VOCs included:

Lab. Project/SDG No.: JC19023 Sample matrix: Soil

No. of Samples: 4

Trip blank No.: -

Field blank No.: -

Equipment blank No.: JC19023-3

Field duplicate No.: -

<input checked="" type="checkbox"/> Data Completeness	<input checked="" type="checkbox"/> Laboratory Control Spikes
<input checked="" type="checkbox"/> Holding Times	<input checked="" type="checkbox"/> Field Duplicates
<input type="checkbox"/> N/A GC/MS Tuning	<input checked="" type="checkbox"/> Calibrations
<input type="checkbox"/> N/A Internal Standard Performance	<input checked="" type="checkbox"/> Compound Identifications
<input checked="" type="checkbox"/> Blanks	<input checked="" type="checkbox"/> Compound Quantitation
<input checked="" type="checkbox"/> Surrogate Recoveries	<input checked="" type="checkbox"/> Quantitation Limits
<input checked="" type="checkbox"/> Matrix Spike/Matrix Spike Duplicate	

Overall Comments: Low molecular weight alcohols by SW-846 8015C

Definition of Qualifiers:

J- Estimated results
 U- Compound not detected
 R- Rejected data
 UJ- Estimated nondetect

Reviewer: Rafael Infante

Date: May 16, 2016

DATA REVIEW WORKSHEETS

DATA COMPLETENESS

MISSING INFORMATION

DATE LAB. CONTACTED

DATE RECEIVED

A hand-drawn graph on lined paper showing a linear relationship between two variables. The x-axis is labeled "Time (hours)" and the y-axis is labeled "Distance (miles)". A line starts at the origin (0,0) and passes through points (1, 50), (2, 100), (3, 150), (4, 200), (5, 250), (6, 300), (7, 350), (8, 400), (9, 450), and (10, 500). The line is labeled "Constant Speed".

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

HOLDING TIMES

The objective of this parameter is to ascertain the validity of the results based on the holding time of the sample from time of collection to the time of analysis.

Complete table for all samples and note the analysis and/or preservation not within criteria

SAMPLE ID	DATE SAMPLED	DATE ANALYZED	pH	ACTION
All samples analyzed within the recommended method holding time. All samples properly preserved except for the followings: Samples JC19023-1; JC19023-2; and JC19023-4: were received in a bulk container but was not preserved within 48 hours of sampling. Analysis performed at Accutest Laboratories, Orlando FL. Results qualified as (UJ) in affected samples.				

Criteria

Aqueous samples – 14 days from sample collection for preserved samples (pH ≤ 2, 4°C), no air bubbles.

Aqueous samples – 7 days from sample collection for unpreserved samples, 4°C, no air bubbles.

Soil samples- 7 days from sample collection.

Cooler temperature (Criteria: 4 ± 2 °C): 5.4°C

Actions

If the VOCs vial(s) have air bubbles, estimate positive results (J) and reject nondetects (R).

If the % solids of soil samples is 10-50%, estimate positive results (J) and nondetects (UJ)

If the % solid of soil samples is < 10%, estimate positive results (J) and reject nondetects (R).

If holding times are exceeded but < 14 days beyond criteria, estimate positive results (J) and nondetects (UJ).

If holding times are exceeded but < 28 days beyond criteria, estimate positive results (J) and reject nondetects (R).

If holding times are grossly exceeded (> 28 days beyond criteria), reject all results (R).

If samples were not iced or if the ice were melted (> 10°C), estimate positive results (J) and nondetects (UJ).

DATA REVIEW WORKSHEETS

All criteria were met ___N/A___
Criteria were not met see below _____

GC/MS TUNING

The assessment of the tuning results is to determine if the sample instrumentation is within the standard tuning QC limits

___N/A___ The BFB performance results were reviewed and found to be within the specified criteria.

___N/A___ BFB tuning was performed for every 12 hours of sample analysis.

If no, use professional judgment to determine whether the associated data should be accepted, qualified or rejected.

List the samples affected: _____

If mass calibration is in error, all associated data are rejected.

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

CALIBRATION VERIFICATION

Compliance requirements for satisfactory instrument calibration are established to ensure that the instrument is capable of producing and maintaining acceptable quantitative data.

Date of initial calibration: 05/02/16
 Dates of continuing calibration: 05/02/16 (initial); 05/03/16
 Dates of final calibration verification: 05/02/16; 05/03/16
 Instrument ID number: VOA5
 Matrix/Level: Aqueous/low

DATE	LAB FILE ID#	CRITERIA OUT RFs, %RSD, %D, r	COMPOUND	SAMPLES AFFECTED

Note: Initial, continuing, and final calibration verifications meets method specific criteria.

Criteria

All RFs must be > 0.05 regardless of method requirements for SPCC.

All %RSD must be $\leq 15\%$ regardless of method requirements for CCC.

All %Ds must be $\leq 20\%$ regardless of method requirements for CCC.

It should be noted that Region 2 SOP HW-24 does not specify criterion for the curve correlation coefficient (r). A limit for r of ≥ 0.995 has therefore been utilized as professional judgment.

Actions

If any compound has an initial RF or a continuing RF of < 0.05 , estimate positive results (J) and reject nondetects (R), regardless of method requirements.

If any compound has a %RSD $> 15\%$, estimate positive results (J) and use professional judgment to qualify nondetects.

If any compound has a %RSD $> 90\%$, estimate positive results (J) and reject nondetects (R).

If any compound has a % D $> 20\%$, estimate positive results (J) and reject nondetects (R).

If any compound has a % D $> 20\%$, estimate positive results (J) and nondetects (UJ).

If any compound has a % D $> 90\%$, estimate positive results (J) and reject nondetects (R).

If any compound has r < 0.995 , estimate positive results and nondetects.

A separate worksheet should be filled for each initial curve

DATA REVIEW WORKSHEETS

All criteria were met X
Criteria were not met
and/or see below _____

V A. BLANK ANALYSIS RESULTS (Sections 1 & 2)

The assessment of the blank analysis results is to determine the existence and magnitude of contamination problems. The criteria for evaluation of blanks apply only to blanks associated with the samples, including trip, equipment, and laboratory blanks. If problems with any blanks exist, all data associated with the case must be carefully evaluated to determine whether or not there is an inherent variability in the data for the case, or if the problem is an isolated occurrence not affecting other data.

List the contamination in the blanks below. High and low levels blanks must be treated separately.

Laboratory blanks

DATE ANALYZED	LAB ID	LEVEL/MATRIX	COMPOUND	CONCENTRATION UNITS
	All_method_blank_meeth_method_specific_criteria			

Field/Equipment/Trip blank

[illegible]

DATA REVIEW WORKSHEETS

All criteria were met X
Criteria were not met
and/or see below _____

V B. BLANK ANALYSIS RESULTS (Section 3)

Blank Actions

Action Levels (ALs) should be based upon the highest concentration of contaminant determined in any blank. Do not qualify any blank with another blank. The ALs for samples which have been diluted should be corrected for the sample dilution factor and/or % moisture, where applicable. No positive sample results should be reported unless the concentration of the compound in the samples exceeds the ALs:

ALs = 10x the amount of common contaminants (methylene chloride, acetone, 2-butanone, and toluene)

ALs = 5x for any other compounds

Specific actions are as follows:

If the concentration is < sample quantitation limit (SQL) and \leq AL, report the compound as not detected (U) at the SQL.

If the concentration is \geq SQL but \leq AL, report the compound as not detected (U) at the reported concentration.

If the concentration is \geq SQL and $>$ AL, report the concentration unqualified.

Notes:

High and low level blanks must be treated separately

Compounds qualified "U" for blank contamination are still considered "hits" when qualifying for calibration criteria.

[illegible]

DATA REVIEW WORKSHEETS

All criteria were met ☒ X
 Criteria were not met
 and/or see below _____

SURROGATE SPIKE RECOVERIES

Laboratory performance of individual samples is established by evaluation of surrogate spike recoveries. All samples are spiked with surrogate compounds prior to sample analysis. The accuracy of the analysis is measured by the surrogate percent recovery. Since the effects of the sample matrix are frequently outside the control of the laboratory and may present relatively unique problems, the validation of data is frequently subjective and demands analytical experience and professional judgment.

List the percent recoveries (%Rs) which do not meet the criteria for surrogate recovery.

Matrix: solid/aqueous

SAMPLE ID	SURROGATE COMPOUND			ACTION
	Hexanol	DBFM	TOL-d8	

____ All surrogate recoveries within laboratory control limits. _____

QC Limits* (Aqueous)

____ LL to UL ____ 73 to 123 ____ to ____ to ____ to ____

QC Limits* (Solid-Low)

____ LL to UL ____ 69 to 121 ____ to ____ to ____ to ____

QC Limits* (Solid-Med)

____ LL to UL ____ to ____ to ____ to ____ to ____

1,2-DCA = 1,2-Dichloromethane-d4

TOL-d8 = Toluene-d8

DBFM = Dibromofluoromethane

BFB = Bromofluorobenzene

* QC limits are laboratory in-house performance criteria, LL = lower limit, UL = upper limit.

* If QC limits are not available, use limits of 80 – 120 % for aqueous and 70 – 130 % for solid samples.

Actions:

QUALITY	%R < 10%	%R = 10% - LL	%R > UL
Positive results	J	J	J
Nondetects results	R	UJ	Accept

Surrogate action should be applied:

If one or more surrogate in the VOC fraction is out of specification, but has a recovery of > 10%.

If any one surrogate in a fraction shows < 10 % recovery.

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

VII. A MATRIX SPIKE/MATRIX SPIKE DUPLICATE (MS/MSD)

This data is generated to determine long term precision and accuracy in the analytical method for various matrices. This data alone cannot be used to evaluate the precision and accuracy of individual samples. If any % R in the MS or MSD falls outside the designated range, the reviewer should determine if there are matrix effects, i.e. LCS data are within the QC limits but MS/MSD data are outside QC limit.

1. MS/MSD Recoveries and Precision Criteria

The laboratory should use one MS and a duplicate analysis of an unspiked field sample if target analytes are expected in the sample. If target analytes are not expected, MS/MSD should be analyzed.

List the %Rs, RPD of the compounds which do not meet the criteria.

Sample ID: JC19214-2AMS/-1MSD Matrix/Level: Aqueous
 Sample ID: JC18972-2AMS/-1MSD Matrix/Level: Soil

MS OR MSD COMPOUND % R RPD QC LIMITS ACTION

MS/MSD%_recoveries_and_RPD_within_laboratory_control_limits_except_for_the_followings:_

The QC reported here applies to the following samples:
 JC19023-1, JC19023-2, JC19023-4

Method: SW846 8015C MOD

Compound	JC18972-2A mg/kg	Q	Spike mg/kg	MS mg/kg	MS %	Spike mg/kg	MSD mg/kg	MSD %	RPD	Limits Rec/RPD
Ethanol	ND		244	284	116	244	305	125*	7	80-117/13
Isobutyl Alcohol	ND		244	278	114	244	294	121*	6	72-117/14
Isopropyl Alcohol	ND		244	294	121*	244	314	129*	7	75-116/15
n-Propyl Alcohol	ND		244	289	118*	244	310	127*	7	78-116/13
n-Butyl Alcohol	ND		244	291	119*	244	310	127*	6	74-115/13
Methanol	ND		244	282	116	244	307	126*	8	77-116/13

Surrogate Recoveries	MS	MSD	JC18972-2A	Limits
Hexanol	115%	120%	107%	69-121%

(a) Sample was received in a bulk container but was not preserved within 48 hours of sampling.

- * QC limits are laboratory in-house performance criteria, LL = lower limit, UL = upper limit.
- * If QC limits are not available, use limits of 70 – 130 %.

Actions:

QUALITY	%R < LL	%R > UL
Positive results	J	J
Nondetects results	R	Accept

All criteria were met X
Criteria were not met
and/or see below

Note: No action taken, MS/MSD applies only the unspiked sample. Unspiked sample from another data package, used for QC purposes only

VII. B MATRIX SPIKE/MATRIX SPIKE DUPLICATE

It should be noted that Region 2 SOP HW-24 does not specify a MS/MSD criteria for the unspiked compounds in the sample. A %RSD of < 50% has therefore been utilized as professional judgment.

List the %RSD of the compounds which do not meet the criteria.

COMPOUND	SAMPLE CONC.	MS CONC.	MSD CONC.	% RSD	ACTION
----------	-----------------	----------	-----------	-------	--------

[illegible]

* If the % RSD > 50, qualify the positive result in the unspiked samples as estimated (J).
* If the % RSD is not calculated (NC) due to nondetected value, use professional judgment to qualify the data.

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

VIII. LABORATORY CONTROL SAMPLE (LCS) ANALYSIS

This data is generated to determine accuracy of the analytical method for various matrices.

1. LCS Recoveries Criteria

Where LCS spiked with the same analyte at the same concentrations as the MS/MSD? **Yes**
 or No. If no make note in data review memo.

List the %R of compounds which do not meet the criteria

LCS ID	COMPOUND	% R	QC LIMIT
Recoveries within laboratory control limits. _____			

- * QC limits are laboratory in-house performance criteria, LL = lower limit, UL = upper limit.
- * If QC limits are not available, use limits of 70 – 130 %.

Actions:

QUALITY	%R < LL	%R > UL
Positive results	J	J
Nondetects results	R	Accept

All analytes in the associated sample results are qualified for the following criteria.

If 25 % of the LCS recoveries were < LL (or 70 %), qualify all positive results (j) and reject nondetects (R).

If two or more LCS were below 10 %, qualify all positive results as (J) and reject nondetects (R).

2. Frequency Criteria:

Where LCS analyzed at the required frequency and for each matrix? Yes or No.

If no, the data may be affected. Use professional judgment to determine the severity of the effect and qualify data accordingly. Discuss any actions below and list the samples affected.

DATA REVIEW WORKSHEETS

All criteria were met N/A
 Criteria were not met
 and/or see below

IX. FIELD/LABORATORY DUPLICATE PRECISION

Sample IDs: - Matrix: -

Field/laboratory duplicate samples may be taken and analyzed as an indication of overall precision. These analyses measure both field and lab precision; therefore, the results may have more variability than laboratory duplicates which only laboratory performance. It is also expected that soil duplicate results will have a greater variance than water matrices due to difficulties associated with collecting identical field duplicate samples.

The project QAPP should be reviewed for project-specific information.

Suggested criteria: RPD \pm 30% for aqueous samples, RPD \pm 50 % for solid samples. If both samples and duplicate are <5 SQL, the RPD criteria is doubled.

COMPOUND	SQL	SAMPLE CONC.	DUPLICATE CONC.	RPD	ACTION
No field/laboratory duplicate analyzed with this data package. MS/MSD % recoveries RPD used to assess precision. RPD within laboratory and generally acceptable control limits.					

Actions:

Qualify as estimated positive results (J) and nondetects (UJ) for the compound that exceeded the above criteria. For organics, only the sample and duplicate will be qualified.

If an RPD cannot be calculated because one or both of the sample results is not detected, the following actions apply:

If one sample result is not detected and the other is greater than 5x the SQL qualify (J/UJ).

If one sample value is not detected and the other is greater than 5x the SQL and the SQLs for the sample and duplicate are significantly different, use professional judgment to determine if qualification is appropriate.

If one sample value is not detected and the other is less than 5x, use professional judgment to determine if qualification is appropriate.

If both sample and duplicate results are not detected, no action is needed.

All criteria were met N/A

Criteria were not met

and/or see below

X. INTERNAL STANDARD PERFORMANCE

The assessment of the internal standard (IS) parameter is used to assist the data reviewer in determining the condition of the analytical instrumentation.

List the internal standard area of samples which do not meet the criteria.

- * Area of +100% or -50% of the IS area in the associated calibration standard.
* Retention time (RT) within 30 seconds of the IS area in the associated calibration standard.

DATE	SAMPLE ID	IS OUT	IS AREA	ACCEPTABLE RANGE	ACTION
------	-----------	--------	---------	------------------	--------

This image shows a single sheet of white paper with horizontal blue or grey ruling lines. The lines are evenly spaced and run across the width of the page. There is no handwriting or other markings on the paper.

Actions:

1. IS actions should be applied to the compound quantitated with the out-of-control ISs

QUALITY	IS AREA < -25%	IS AREA = -25 % TO - 50%	IS AREA > + 100%
Positive results	J	J	J
Nondetected results	R	UJ	ACCEPT

2. If a IS retention time varies more than 30 seconds, the chromatographic profile for that sample must be examined to determine if any false positive or negative exists. For shifts of a large magnitude, the reviewer may consider partial or total rejection of the data for the sample fraction.

DATA REVIEW WORKSHEETS

All criteria were met X
Criteria were not met
and/or see below

XII. SAMPLE QUANTITATION

The sample quantitation evaluation is to verify laboratory quantitation results. In the space below, please show a minimum of one sample calculation:

JC19023-1

Hexanol

RF = 5304

$$[] = (587600)/(5304)$$

$$= 110.8 \text{ ppm OK}$$

All criteria were met X

Criteria were not met

and/or see below _____

XII. QUANTITATION LIMITS

A. Dilution performed

[illegible]

B. Percent Solids

List samples which have $\leq 50\%$ solids

[Blank lined paper area]

Actions:

If the % solids of a soil sample is 10-50%, estimate positive results (J) and nondetects (UJ)

If the % solids of a soil sample is $< 10\%$, estimate positive results (J) and reject nondetects (R)

EXECUTIVE NARRATIVE

SDG No: **JC19214** Laboratory: **Accutest, New Jersey**
Analysis: **SW846-8270D** Number of Samples: **4**
Location: **BMSMC, Building 5 Area**
Humacao, PR

SUMMARY: Four (4) samples were analyzed for the ABN TCL list following method SW846-8270D; Naphthalene and 1,4-Dioxane were also analyzed by SW846-8270D using the selective ion monitoring (SIM) technique. The sample results were assessed according to USEPA data validation guidance documents in the following order of precedence: EPA Hazardous Waste Support Section, SOP HW-35A, July 2015 –Revision 0. *Semivolatile Data Validation*. The QC criteria and data validation actions listed on the data review worksheets are from the primary guidance document, unless otherwise noted.

Results are valid and can be used for decision making purposes.

Critical issues: **None**
Major: **None**
Minor: **None**

Critical findings: **None**
Major findings: **None**
Minor findings: 1. Initial and continuing calibration verifications meet the required criteria. Analytes not meeting the method % difference criteria meet the guidance document performance criteria for continuing calibration verification of ± 25 or 40 %, no action taken. No closing calibration verification included in data package. No action taken, professional judgment.
2. Analytes not meeting the continuing calibration verification criteria of the guidance document were qualified UJ in sample JC19023-1.
3. MS/MSD recoveries for 1,4-Dioxane in sample JC19023-1 outside the laboratory lower recovery limit. Results for 1,4-Dioxane qualified as estimated (J or UJ) in affected samples.

COMMENTS: Results are valid and can be used for decision making purposes.

Reviewers Name: **Rafael Infante**
Chemist License 1888

Signature:



Date: **May 17, 2016**

SAMPLE ORGANIC DATA SAMPLE SUMMARY

Sample ID: JC19023-1
Sample location: BMSMC Building 5 Area
Sampling date: 4/22/2016
Matrix: Soil

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
2-Chlorophenol	88	ug/kg	1	-	U	Yes
4-Chloro-3-methyl phenol	220	ug/kg	1	-	U	Yes
2,4-Dichlorophenol	220	ug/kg	1	-	U	Yes
2,4-Dimethylphenol	220	ug/kg	1	-	U	Yes
2,4-Dinitrophenol	220	ug/kg	1	-	UJ	Yes
4,6-Dinitro-o-cresol	220	ug/kg	1	-	U	Yes
2-Methylphenol	88	ug/kg	1	-	U	Yes
3&4-Methylphenol	88	ug/kg	1	-	U	Yes
2-Nitrophenol	220	ug/kg	1	-	U	Yes
4-Nitrophenol	450	ug/kg	1	-	U	Yes
Pentachlorophenol	220	ug/kg	1	-	UJ	Yes
Phenol	88	ug/kg	1	-	U	Yes
2,3,4,6-Tetrachlorophenol	220	ug/kg	1	-	U	Yes
2,4,5-Trichlorophenol	220	ug/kg	1	-	U	Yes
2,4,6-Trichlorophenol	220	ug/kg	1	-	U	Yes
Acenaphthene	45	ug/kg	1	-	U	Yes
Acenaphthylene	45	ug/kg	1	-	U	Yes
Acetophenone	220	ug/kg	1	-	U	Yes
Anthracene	45	ug/kg	1	-	U	Yes
Atrazine	88	ug/kg	1	-	UJ	Yes
Benzo(a)anthracene	45	ug/kg	1	-	U	Yes
Benzo(a)pyrene	45	ug/kg	1	-	U	Yes
Benzo(b)fluoranthene	45	ug/kg	1	-	U	Yes
Benzo(g,h,i)perylene	45	ug/kg	1	-	U	Yes
Benzo(k)fluoranthene	45	ug/kg	1	-	U	Yes
4-Bromophenyl phenyl ether	88	ug/kg	1	-	U	Yes
Butyl benzyl phthalate	88	ug/kg	1	-	U	Yes
1,1'-Biphenyl	88	ug/kg	1	-	U	Yes
Benzaldehyde	103	ug/kg	1	J	UJ	Yes
2-Chloronaphthalene	88	ug/kg	1	-	U	Yes
4-Chloroaniline	220	ug/kg	1	-	U	Yes
Carbazole	88	ug/kg	1	-	U	Yes
Caprolactam	88	ug/kg	1	-	U	Yes
Chrysene	45	ug/kg	1	-	U	Yes
bis(2-Chloroethoxy)methane	88	ug/kg	1	-	U	Yes
bis(2-Chloroethyl)ether	88	ug/kg	1	-	U	Yes

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
bis(2-Chloroisopropyl)ether	88	ug/kg	1	-	U	Yes
4-Chlorophenyl phenyl ether	88	ug/kg	1	-	U	Yes
2,4-Dinitrotoluene	45	ug/kg	1	-	U	Yes
2,6-Dinitrotoluene	45	ug/kg	1	-	U	Yes
3,3'-Dichlorobenzidine	88	ug/kg	1	-	U	Yes
Dibenzo(a,h)anthracene	45	ug/kg	1	-	U	Yes
Dibenzofuran	88	ug/kg	1	-	U	Yes
Di-n-butyl phthalate	88	ug/kg	1	-	U	Yes
Di-n-octyl phthalate	88	ug/kg	1	-	U	Yes
Diethyl phthalate	88	ug/kg	1	-	U	Yes
Dimethyl phthalate	88	ug/kg	1	-	U	Yes
bis(2-Ethylhexyl)phthalate	88	ug/kg	1	-	U	Yes
Fluoranthene	24.3	ug/kg	1	J	UJ	Yes
Fluorene	45	ug/kg	1	-	U	Yes
Hexachlorobenzene	88	ug/kg	1	-	U	Yes
Hexachlorobutadiene	45	ug/kg	1	-	U	Yes
Hexachlorocyclopentadiene	450	ug/kg	1	-	U	Yes
Hexachloroethane	220	ug/kg	1	-	U	Yes
Indeno(1,2,3-cd)pyrene	45	ug/kg	1	-	U	Yes
Isophorone	88	ug/kg	1	-	U	Yes
1-Methylnaphthalene	88	ug/kg	1	-	U	Yes
2-Methylnaphthalene	88	ug/kg	1	-	U	Yes
2-Nitroaniline	220	ug/kg	1	-	U	Yes
3-Nitroaniline	220	ug/kg	1	-	U	Yes
4-Nitroaniline	220	ug/kg	1	-	U	Yes
Nitrobenzene	88	ug/kg	1	-	U	Yes
N-Nitroso-di-n-propylamine	88	ug/kg	1	-	U	Yes
Nitrosodiphenylamine	220	ug/kg	1	-	U	Yes
Phenanthrene	33.9	ug/kg	1	J	UJ	Yes
Pyrene	19.1	ug/kg	1	J	UJ	Yes
1,2,4,5-Tetrachlorobenzene	220	ug/kg	1	-	U	Yes

METHOD: 8270D (SIM)

Naphthalene	30.0	ug/kg	1	-	-	Yes
1,4-Dioxane	7.80	ug/kg	1	-	J	Yes

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
--------------	--------	-------	-----------------	----------	------------	------------

Sample ID: JC19023-2

Sample location: BMSMC Building 5 Area

Sampling date: 4/22/2016

Matrix: Soil

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
2-Chlorophenol	89	ug/Kg	1	-	U	Yes
4-Chloro-3-methyl phenol	220	ug/Kg	1	-	U	Yes
2,4-Dichlorophenol	220	ug/Kg	1	-	U	Yes
2,4-Dimethylphenol	220	ug/Kg	1	-	U	Yes
2,4-Dinitrophenol	220	ug/Kg	1	-	U	Yes
4,6-Dinitro-o-cresol	220	ug/Kg	1	-	U	Yes
2-Methylphenol	89	ug/Kg	1	-	U	Yes
3&4-Methylphenol	89	ug/Kg	1	-	U	Yes
2-Nitrophenol	220	ug/Kg	1	-	U	Yes
4-Nitrophenol	450	ug/Kg	1	-	U	Yes
Pentachlorophenol	220	ug/Kg	1	-	U	Yes
Phenol	89	ug/Kg	1	-	U	Yes
2,3,4,6-Tetrachlorophenol	220	ug/Kg	1	-	U	Yes
2,4,5-Trichlorophenol	220	ug/Kg	1	-	U	Yes
2,4,6-Trichlorophenol	220	ug/Kg	1	-	U	Yes
Acenaphthene	45	ug/Kg	1	-	U	Yes
Acenaphthylene	45	ug/Kg	1	-	U	Yes
Acetophenone	220	ug/Kg	1	-	U	Yes
Anthracene	45	ug/Kg	1	-	U	Yes
Atrazine	89	ug/Kg	1	-	U	Yes
Benzo(a)anthracene	45	ug/Kg	1	-	U	Yes
Benzo(a)pyrene	45	ug/Kg	1	-	U	Yes
Benzo(b)fluoranthene	45	ug/Kg	1	-	U	Yes
Benzo(g,h,i)perylene	45	ug/Kg	1	-	U	Yes
Benzo(k)fluoranthene	45	ug/Kg	1	-	U	Yes
4-Bromophenyl phenyl ether	89	ug/Kg	1	-	U	Yes
Butyl benzyl phthalate	89	ug/Kg	1	-	U	Yes
1,1'-Biphenyl	89	ug/Kg	1	-	U	Yes
Benzaldehyde	220	ug/Kg	1	-	U	Yes
2-Chloronaphthalene	89	ug/Kg	1	-	U	Yes
4-Chloroaniline	220	ug/Kg	1	-	U	Yes
Carbazole	89	ug/Kg	1	-	U	Yes
Caprolactam	89	ug/Kg	1	-	U	Yes
Chrysene	45	ug/Kg	1	-	U	Yes
bis(2-Chloroethoxy)methane	89	ug/Kg	1	-	U	Yes
bis(2-Chloroethyl)ether	89	ug/Kg	1	-	U	Yes

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
bis(2-Chloroisopropyl)ether	89	ug/Kg	1	-	U	Yes
4-Chlorophenyl phenyl ether	89	ug/Kg	1	-	U	Yes
2,4-Dinitrotoluene	45	ug/Kg	1	-	U	Yes
2,6-Dinitrotoluene	45	ug/Kg	1	-	U	Yes
3,3'-Dichlorobenzidine	89	ug/Kg	1	-	U	Yes
Dibenzo(a,h)anthracene	45	ug/Kg	1	-	U	Yes
Dibenzofuran	89	ug/Kg	1	-	U	Yes
Di-n-butyl phthalate	89	ug/Kg	1	-	U	Yes
Di-n-octyl phthalate	89	ug/Kg	1	-	U	Yes
Diethyl phthalate	89	ug/Kg	1	-	U	Yes
Dimethyl phthalate	89	ug/Kg	1	-	U	Yes
bis(2-Ethylhexyl)phthalate	89	ug/Kg	1	-	U	Yes
Fluoranthene	45	ug/Kg	1	-	U	Yes
Fluorene	45	ug/Kg	1	-	U	Yes
Hexachlorobenzene	89	ug/Kg	1	-	U	Yes
Hexachlorobutadiene	45	ug/Kg	1	-	U	Yes
Hexachlorocyclopentadiene	450	ug/Kg	1	-	U	Yes
Hexachloroethane	220	ug/Kg	1	-	U	Yes
Indeno(1,2,3-cd)pyrene	45	ug/Kg	1	-	U	Yes
Isophorone	89	ug/Kg	1	-	U	Yes
1-Methylnaphthalene	89	ug/Kg	1	-	U	Yes
2-Methylnaphthalene	89	ug/Kg	1	-	U	Yes
2-Nitroaniline	220	ug/Kg	1	-	U	Yes
3-Nitroaniline	220	ug/Kg	1	-	U	Yes
4-Nitroaniline	220	ug/Kg	1	-	U	Yes
Nitrobenzene	89	ug/Kg	1	-	U	Yes
N-Nitroso-di-n-propylamine	89	ug/Kg	1	-	U	Yes
Nitrosodiphenylamine	220	ug/Kg	1	-	U	Yes
Phenanthrene	45	ug/Kg	1	-	U	Yes
Pyrene	45	ug/Kg	1	-	U	Yes
1,2,4,5-Tetrachlorobenzene	220	ug/Kg	1	-	U	Yes

METHOD: 8270D (SIM)

Naphthalene	4.5	ug/Kg	1	-	U	Yes
1,4-Dioxane	4.5	ug/Kg	1	-	UJ	Yes

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
--------------	--------	-------	-----------------	----------	------------	------------

Sample ID: JC19023-3

Sample location: BMSMC Building 5 Area

Sampling date: 4/22/2016

Matrix: AQ - Equipment Blank

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
2-Chlorophenol	5.0	ug/l	1	-	U	Yes
4-Chloro-3-methyl phenol	5.0	ug/l	1	-	U	Yes
2,4-Dichlorophenol	2.0	ug/l	1	-	U	Yes
2,4-Dimethylphenol	5.0	ug/l	1	-	U	Yes
2,4-Dinitrophenol	10	ug/l	1	-	U	Yes
4,6-Dinitro-o-cresol	5.0	ug/l	1	-	U	Yes
2-Methylphenol	2.0	ug/l	1	-	U	Yes
3&4-Methylphenol	2.0	ug/l	1	-	U	Yes
2-Nitrophenol	5.0	ug/l	1	-	U	Yes
4-Nitrophenol	10	ug/l	1	-	U	Yes
Pentachlorophenol	5.0	ug/l	1	-	U	Yes
Phenol	2.0	ug/l	1	-	U	Yes
2,3,4,6-Tetrachlorophenol	5.0	ug/l	1	-	U	Yes
2,4,5-Trichlorophenol	5.0	ug/l	1	-	U	Yes
2,4,6-Trichlorophenol	5.0	ug/l	1	-	U	Yes
Acenaphthene	1.0	ug/l	1	-	U	Yes
Acenaphthylene	1.0	ug/l	1	-	U	Yes
Acetophenone	2.0	ug/l	1	-	U	Yes
Anthracene	1.0	ug/l	1	-	U	Yes
Atrazine	2.0	ug/l	1	-	U	Yes
Benzaldehyde	5.0	ug/l	1	-	U	Yes
Benzo(a)anthracene	1.0	ug/l	1	-	U	Yes
Benzo(a)pyrene	1.0	ug/l	1	-	U	Yes
Benzo(b)fluoranthene	1.0	ug/l	1	-	U	Yes
Benzo(g,h,i)perylene	1.0	ug/l	1	-	U	Yes
Benzo(k)fluoranthene	1.0	ug/l	1	-	U	Yes
4-Bromophenyl phenyl ether	5.0	ug/l	1	-	U	Yes
Butyl benzyl phthalate	5.0	ug/l	1	-	U	Yes
1,1'-Biphenyl	5.0	ug/l	1	-	U	Yes
2-Chloronaphthalene	5.0	ug/l	1	-	U	Yes
4-Chloroaniline	220.0	ug/l	1	-	U	Yes
Carbazole	1.0	ug/l	1	-	U	Yes
Caprolactam	2.0	ug/l	1	-	U	Yes
Chrysene	1.0	ug/l	1	-	U	Yes
bis(2-Chloroethoxy)methane	2.0	ug/l	1	-	U	Yes
bis(2-Chloroethyl)ether	2.0	ug/l	1	-	U	Yes

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
bis(2-Chloroisopropyl)ether	2.0	ug/l	1	-	U	Yes
4-Chlorophenyl phenyl ether	2.0	ug/l	1	-	U	Yes
2,4-Dinitrotoluene	1.0	ug/l	1	-	U	Yes
2,6-Dinitrotoluene	1.0	ug/l	1	-	U	Yes
3,3'-Dichlorobenzidine	2.0	ug/l	1	-	U	Yes
Dibenzo(a,h)anthracene	1.0	ug/l	1	-	U	Yes
Dibenzofuran	5.0	ug/l	1	-	U	Yes
Di-n-butyl phthalate	1.8	ug/l	1	JB	UJ	Yes
Di-n-octyl phthalate	2.0	ug/l	1	-	U	Yes
Diethyl phthalate	2.0	ug/l	1	-	U	Yes
Dimethyl phthalate	2.0	ug/l	1	-	U	Yes
bis(2-Ethylhexyl)phthalate	2.0	ug/l	1	-	U	Yes
Fluoranthene	1.0	ug/l	1	-	U	Yes
Fluorene	1.0	ug/l	1	-	U	Yes
Hexachlorobenzene	1.0	ug/l	1	-	U	Yes
Hexachlorobutadiene	1.0	ug/l	1	-	U	Yes
Hexachlorocyclopentadiene	10	ug/l	1	-	U	Yes
Hexachloroethane	2.0	ug/l	1	-	U	Yes
Indeno(1,2,3-cd)pyrene	1.0	ug/l	1	-	U	Yes
Isophorone	2.0	ug/l	1	-	U	Yes
1-Methylnaphthalene	1.0	ug/l	1	-	U	Yes
2-Methylnaphthalene	1.0	ug/l	1	-	U	Yes
2-Nitroaniline	5.0	ug/l	1	-	U	Yes
3-Nitroaniline	5.0	ug/l	1	-	U	Yes
4-Nitroaniline	5.0	ug/l	1	-	U	Yes
Nitrobenzene	2.0	ug/l	1	-	U	Yes
N-Nitroso-di-n-propylamine	2.0	ug/l	1	-	U	Yes
Nitrosodiphenylamine	5.0	ug/l	1	-	U	Yes
Phenanthrene	1.0	ug/l	1	-	U	Yes
Pyrene	1.0	ug/l	1	-	U	Yes
1,2,4,5-Tetrachlorobenzene	2.0	ug/l	1	-	U	Yes

METHOD: 8270D (SIM)

Naphthalene	0.10	ug/l	1	-	U	Yes
1,4-Dioxane	0.10	ug/l	1	-	U	Yes

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
--------------	--------	-------	-----------------	----------	------------	------------

Sample ID: JC19023-4

Sample location: BMSMC Building 5 Area

Sampling date: 4/25/2016

Matrix: Soil

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
2-Chlorophenol	85	ug/kg	1	-	U	Yes
4-Chloro-3-methyl phenol	210	ug/kg	1	-	U	Yes
2,4-Dichlorophenol	210	ug/kg	1	-	U	Yes
2,4-Dimethylphenol	210	ug/kg	1	-	U	Yes
2,4-Dinitrophenol	210	ug/kg	1	-	U	Yes
4,6-Dinitro-o-cresol	210	ug/kg	1	-	U	Yes
2-Methylphenol	85	ug/kg	1	-	U	Yes
3&4-Methylphenol	85	ug/kg	1	-	U	Yes
2-Nitrophenol	210	ug/kg	1	-	U	Yes
4-Nitrophenol	420	ug/kg	1	-	U	Yes
Pentachlorophenol	210	ug/kg	1	-	U	Yes
Phenol	85	ug/kg	1	-	U	Yes
2,3,4,6-Tetrachlorophenol	210	ug/kg	1	-	U	Yes
2,4,5-Trichlorophenol	210	ug/kg	1	-	U	Yes
2,4,6-Trichlorophenol	210	ug/kg	1	-	U	Yes
Acenaphthene	42	ug/kg	1	-	U	Yes
Acenaphthylene	42	ug/kg	1	-	U	Yes
Acetophenone	210	ug/kg	1	-	U	Yes
Anthracene	72.1	ug/kg	1	-	-	Yes
Atrazine	85	ug/kg	1	-	U	Yes
Benzo(a)anthracene	42	ug/kg	1	-	U	Yes
Benzo(a)pyrene	42	ug/kg	1	-	U	Yes
Benzo(b)fluoranthene	42	ug/kg	1	-	U	Yes
Benzo(g,h,i)perylene	42	ug/kg	1	-	U	Yes
Benzo(k)fluoranthene	85	ug/kg	1	-	U	Yes
4-Bromophenyl phenyl ether	85	ug/kg	1	-	U	Yes
Butyl benzyl phthalate	85	ug/kg	1	-	U	Yes
1,1'-Biphenyl	85	ug/kg	1	-	U	Yes
Benzaldehyde	210	ug/kg	1	-	U	Yes
2-Chloronaphthalene	85	ug/kg	1	-	U	Yes
4-Chloroaniline	210	ug/kg	1	-	U	Yes
Carbazole	85	ug/kg	1	-	U	Yes
Caprolactam	85	ug/kg	1	-	U	Yes
Chrysene	42	ug/kg	1	-	U	Yes
bis(2-Chloroethoxy)methane	85	ug/kg	1	-	U	Yes
bis(2-Chloroethyl)ether	85	ug/kg	1	-	U	Yes

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
bis(2-Chloroisopropyl)ether	85	ug/kg	1	-	U	Yes
4-Chlorophenyl phenyl ether	85	ug/kg	1	-	U	Yes
2,4-Dinitrotoluene	42	ug/kg	1	-	U	Yes
2,6-Dinitrotoluene	42	ug/kg	1	-	U	Yes
3,3'-Dichlorobenzidine	85	ug/kg	1	-	U	Yes
Dibenzo(a,h)anthracene	42	ug/kg	1	-	U	Yes
Dibenzofuran	85	ug/kg	1	-	U	Yes
Di-n-butyl phthalate	85	ug/kg	1	-	U	Yes
Di-n-octyl phthalate	85	ug/kg	1	-	U	Yes
Diethyl phthalate	85	ug/kg	1	-	U	Yes
Dimethyl phthalate	85	ug/kg	1	-	U	Yes
bis(2-Ethylhexyl)phthalate	85.0	ug/kg	1	-	U	Yes
Fluoranthene	42	ug/kg	1	-	U	Yes
Fluorene	42	ug/kg	1	-	U	Yes
Hexachlorobenzene	85	ug/kg	1	-	U	Yes
Hexachlorobutadiene	42	ug/kg	1	-	U	Yes
Hexachlorocyclopentadiene	420	ug/kg	1	-	U	Yes
Hexachloroethane	210	ug/kg	1	-	U	Yes
Indeno(1,2,3-cd)pyrene	85	ug/kg	1	-	U	Yes
Isophorone	85	ug/kg	1	-	U	Yes
2-Methylnaphthalene	85	ug/kg	1	-	U	Yes
2-Nitroaniline	210	ug/kg	1	-	U	Yes
3-Nitroaniline	210	ug/kg	1	-	U	Yes
4-Nitroaniline	210	ug/kg	1	-	U	Yes
Nitrobenzene	85	ug/kg	1	-	U	Yes
N-Nitroso-di-n-propylamine	85	ug/kg	1	-	U	Yes
Nitrosodiphenylamine	210	ug/kg	1	-	U	Yes
Phenanthrene	42	ug/kg	1	-	U	Yes
Pyrene	42	ug/kg	1	-	U	Yes
1,2,4,5-Tetrachlorobenzene	210	ug/kg	1	-	U	Yes

METHOD: 8270D (SIM)

Naphthalene	4.2	ug/L	1	-	-	Yes
1,4-Dioxane	4.2	ug/kg	1	-	UJ	Yes

DATA REVIEW WORKSHEETS

Project Number: JC19023

Date: April 22-25, 2016

Shipping Date: April 25, 2016

EPA Region: 2

REVIEW OF SEMIVOLATILE ORGANIC PACKAGE

The following guidelines for evaluating volatile organics were created to delineate required validation actions. This document will assist the reviewer in using professional judgment to make more informed decision and in better serving the needs of the data users. The sample results were assessed according to USEPA data validation guidance documents in the following order of precedence: EPA Hazardous Waste Support Section, SOP HW-35A, July 2015 –Revision 0. *Semivolatile Data Validation*. The QC criteria and data validation actions listed on the data review worksheets are from the primary guidance document, unless otherwise noted.

The hardcopied (laboratory name) Accutest data package received has been reviewed and the quality control and performance data summarized. The data review for SVOCs included:

Lab. Project/SDG No.: JC19023 Sample matrix: Groundwater

No. of Samples: 4 Full_scan/4_SIM

Trip blank No.: -

Field blank No.: -

Equipment blank No.: JC19023-3

Field duplicate No.: -

<input checked="" type="checkbox"/> Data Completeness	<input checked="" type="checkbox"/> Laboratory Control Spikes
<input checked="" type="checkbox"/> Holding Times	<input checked="" type="checkbox"/> Field Duplicates
<input checked="" type="checkbox"/> GC/MS Tuning	<input checked="" type="checkbox"/> Calibrations
<input checked="" type="checkbox"/> Internal Standard Performance	<input checked="" type="checkbox"/> Compound Identifications
<input checked="" type="checkbox"/> Blanks	<input checked="" type="checkbox"/> Compound Quantitation
<input checked="" type="checkbox"/> Surrogate Recoveries	<input checked="" type="checkbox"/> Quantitation Limits
<input checked="" type="checkbox"/> Matrix Spike/Matrix Spike Duplicate	

Overall Comments: ABN TCL list by method SW846-8270D; Naphthalene and 1,4-Dioxane analyzed by method SW846-8270D (SIM)

Definition of Qualifiers:

J- Estimated results
 U- Compound not detected
 R- Rejected data
 UJ- Estimated nondetect

Reviewer: Rafael Infante

Date: May 16, 2016

DATA REVIEW WORKSHEETS

DATA COMPLETENESS

MISSING INFORMATION

DATE LAB. CONTACTED

DATE RECEIVED

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

HOLDING TIMES

The objective of this parameter is to ascertain the validity of the results based on the holding time of the sample from time of collection to the time of analysis.

Complete table for all samples and note the analysis and/or preservation not within criteria

SAMPLE ID	DATE SAMPLED	DATE EXTRACTED/ANALYZED	pH	ACTION
All samples extracted and analyzed within method recommended holding time.				

Cooler temperature (Criteria: 4 ± 2 °C): 4.4°C

Actions

Results will be qualified based on the criteria of the following Table:

Table 1. Holding Time Actions for Semivolatile Analyses

Matrix	Preserved	Criteria	Action	
			Detected Associated Compounds	Non-Detected Associated Compounds
Aqueous	No	≤ 7 days (for extraction) ≤ 40 days (for analysis)	Use professional judgment	
	No	> 7 days (for extraction) > 40 days (for analysis)	J	Use professional judgment
	Yes	≤ 7 days (for extraction) ≤ 40 days (for analysis)	No qualification	
	Yes	> 7 days (for extraction) > 40 days (for analysis)	J	UJ
	Yes/No	Grossly Exceeded	J	UJ or R
Non-Aqueous	No	≤ 14 days (for extraction) ≤ 40 days (for analysis)	Use professional judgment	
	No	> 14 days (for extraction) > 40 days (for analysis)	J	Use professional judgment
	Yes	≤ 14 days (for extraction) ≤ 40 days (for analysis)	No qualification	
	Yes	> 14 days (for extraction) > 40 days (for analysis)	J	UJ
	Yes/No	Grossly Exceeded	J	UJ or R

DATA REVIEW WORKSHEETS

All criteria were met X
Criteria were not met see below

GC/MS TUNING

The assessment of the tuning results is to determine if the sample instrumentation is within the standard tuning QC limits

 X The DFTPP performance results were reviewed and found to be within the specified criteria.

 X DFTPP tuning was performed for every 12 hours of sample analysis.

If no, use professional judgment to determine whether the associated data should be accepted, qualified or rejected.

Notes: These requirements do not apply when samples are analyzed by the Selected Ion Monitoring (SIM) technique.

All mass spectrometer conditions must be identical to those used during the sample analysis. Background subtraction actions resulting in spectral distortion are unacceptable

Notes: No data should be qualified based of DFTPP failure.

The requirement to analyze the instrument performance check solution is optional when analysis of PAHs/pentachlorophenol is to be performed by the SIM technique.

List	the	samples	affected:
_____			_____
_____			_____
_____			_____
_____			_____

Actions:

1. If sample are analyzed without a preceding valid instrument performance check or are analyzed 12 hours after the Instrument Performance Check, qualify all data in those samples as unusable (R).
2. If ion abundance criteria are not met, use professional judgment to determine to what extent the data may be utilized.
3. State in the Data Review Narrative, decisions to use analytical data associated with DFTPP instrument performance checks not meeting the contract requirements.
4. Use professional judgment to determine if associated data should be qualified based on the spectrum of the mass calibration compounds.

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

INITIAL CALIBRATION VERIFICATION

Compliance requirements for satisfactory instrument calibration are established to ensure that the instrument is capable of producing and maintaining acceptable quantitative data.

Date of initial calibration: <u>04/21/2016_(SIM)</u>	<u>04/14/2016_(SIM)</u>
Instrument ID numbers: <u>GCMS3M</u>	<u>GCMS4M</u>
Matrix/Level: <u>Aqueous/low</u>	<u>Aqueous/low</u>
Date of initial calibration: <u>04/19/16; 04/25/16_(Scan)</u>	<u>04/12/16; 04/26/16_(Scan)</u>
Instrument ID numbers: <u>GCMS2P</u>	<u>GCMSM</u>
Matrix/Level: <u>Aqueous/low</u>	<u>Aqueous/low</u>
Date of initial calibration: <u>04/11/16_(Scan)</u>	
Instrument ID numbers: <u>GCMS5P</u>	
Matrix/Level: <u>Aqueous/low</u>	

DATE	LAB ID#	FILE	CRITERIA OUT RFs, %RSD, %D, r	COMPOUND	SAMPLES AFFECTED

Actions:

Qualify the initial calibration analytes listed in Table 2 using the following criteria:

Table 3. Initial Calibration Actions for Semivolatile Analysis

Criteria	Action	
	Detect	Non-detect
Initial Calibration not performed at specified frequency and sequence	Use professional judgment R	Use professional judgment R
Initial Calibration not performed at the specified concentrations	J	UJ
RRF < Minimum RRF in Table 2 for target analyte	Use professional judgment J+ or R	R
RRF ≥ Minimum RRF in Table 2 for target analyte	No qualification	No qualification
%RSD > Maximum %RSD in Table 2 for target analyte	J	Use professional judgment
%RSD ≤ Maximum %RSD in Table 2 for target analyte	No qualification	No qualification

Initial Calibration**Table 2. RRF, %RSD, and %D Acceptance Criteria in Initial Calibration and CCV for Semivolatil Analysis**

Analyte	Minimum RRF	Maximum %RSD	Opening Maximum %D ¹	Opening Maximum %D ¹
1,4-Dioxane	0.010	40.0	± 40.0	± 50.0
Benzaldehyde	0.100	40.0	± 40.0	± 50.0
Phenol	0.080	20.0	± 20.0	± 25.0
Bis(2-chloroethyl)ether	0.100	20.0	± 20.0	± 25.0
2-Chlorophenol	0.200	20.0	± 20.0	± 25.0
2-Methylphenol	0.010	20.0	± 20.0	± 25.0
3-Methylphenol	0.010	20.0	± 20.0	± 25.0
2,2'-Oxybis-(1-chloropropane)	0.010	20.0	± 25.0	± 50.0
Acetophenone	0.060	20.0	± 20.0	± 25.0
4-Methylphenol	0.010	20.0	± 20.0	± 25.0
N-Nitroso-di-n-propylamine	0.080	20.0	± 25.0	± 25.0
Hexachloroethane	0.100	20.0	± 20.0	± 25.0
Nitrobenzene	0.090	20.0	± 20.0	± 25.0
Isophorone	0.100	20.0	± 20.0	± 25.0
2-Nitrophenol	0.060	20.0	± 20.0	± 25.0
2,4-Dimethylphenol	0.050	20.0	± 25.0	± 50.0
Bis(2-chloroethoxy)methane	0.080	20.0	± 20.0	± 25.0
2,4-Dichlorophenol	0.060	20.0	± 20.0	± 25.0
Naphthalene	0.200	20.0	± 20.0	± 25.0
4-Chloroaniline	0.010	40.0	± 40.0	± 50.0
Hexachlorobutadiene	0.040	20.0	± 20.0	± 25.0
Caprolactam	0.010	40.0	± 30.0	± 50.0
4-Chloro-3-methylphenol	0.040	20.0	± 20.0	± 25.0
2-Methylnaphthalene	0.100	20.0	± 20.0	± 25.0
Hexachlorocyclopentadiene	0.010	40.0	± 40.0	± 50.0
2,4,6-Trichlorophenol	0.090	20.0	± 20.0	± 25.0
2,4,5-Trichlorophenol	0.100	20.0	± 20.0	± 25.0
1,1'-Biphenyl	0.200	20.0	± 20.0	± 25.0

DATA REVIEW WORKSHEETS

Analyte	Minimum RRF	Maximum %RSD	Opening Maximum %D ¹	Opening Maximum %D ¹
2-Chloronaphthalene	0.300	20.0	± 20.0	± 25.0
2-Nitroaniline	0.060	20.0	± 25.0	± 25.0
Dimethylphthalate	0.300	20.0	± 25.0	± 25.0
2,6-Dinitrotoluene	0.080	20.0	± 20.0	± 25.0
Acenaphthylene	0.400	20.0	± 20.0	± 25.0
3-Nitroaniline	0.010	20.0	± 25.0	± 50.0
Acenaphthene	0.200	20.0	± 20.0	± 25.0
2,4-Dinitrophenol	0.010	40.0	± 50.0	± 50.0
4-Nitrophenol	0.010	40.0	± 40.0	± 50.0
Dibenzofuran	0.300	20.0	± 20.0	± 25.0
2,4-Dinitrotoluene	0.070	20.0	± 20.0	± 25.0
Diethylphthalate	0.300	20.0	± 20.0	± 25.0
1,2,4,5-Tetrachlorobenzene	0.100	20.0	± 20.0	± 25.0
4-Chlorophenyl-phenylether	0.100	20.0	± 20.0	± 25.0
Fluorene	0.200	20.0	± 20.0	± 25.0
4-Nitroaniline	0.010	40.0	± 40.0	± 50.0
4,6-Dinitro-2-methylphenol	0.010	40.0	± 30.0	± 50.0
4-Bromophenyl-phenyl ether	0.070	20.0	± 20.0	± 25.0
N-Nitrosodiphenylamine	0.100	20.0	± 20.0	± 25.0
Hexachlorobenzene	0.050	20.0	± 20.0	± 25.0
Atrazine	0.010	40.0	± 25.0	± 50.0
Pentachlorophenol	0.010	40.0	± 40.0	± 50.0
Phenanthrene	0.200	20.0	± 20.0	± 25.0
Anthracene	0.200	20.0	± 20.0	± 25.0
Carbazole	0.050	20.0	± 20.0	± 25.0
Di-n-butylphthalate	0.500	20.0	± 20.0	± 25.0
Fluoranthene	0.100	20.0	± 20.0	± 25.0
Pyrene	0.400	20.0	± 25.0	± 50.0
Butylbenzylphthalate	0.100	20.0	± 25.0	± 50.0

DATA REVIEW WORKSHEETS

Analyte	Minimum RRF	Maximum %RSD	Opening Maximum %D ¹	Opening Maximum %D ¹
3,3'-Dichlorobenzidine	0.010	40.0	± 40.0	± 50.0
Benzo(a)anthracene	0.300	20.0	± 20.0	± 25.0
Chrysene	0.200	20.0	± 20.0	± 50.0
Bis(2-ethylhexyl) phthalate	0.200	20.0	± 25.0	± 50.0
Di-n-octylphthalate	0.010	40.0	± 40.0	± 50.0
Benzo(b)fluoranthene	0.010	20.0	± 25.0	± 50.0
Benzo(k)fluoranthene	0.010	20.0	± 25.0	± 50.0
Benzo(a)pyrene	0.010	20.0	± 20.0	± 50.0
Indeno(1,2,3-cd)pyrene	0.010	20.0	± 25.0	± 50.0
Dibenzo(a,h)anthracene	0.010	20.0	± 25.0	± 50.0
Benzo(g,h,i)perylene	0.010	20.0	± 30.0	± 50.0
2,3,4,6-Tetrachlorophenol	0.040	20.0	± 20.0	± 50.0
Naphthalene	0.600	20.0	± 25.0	± 25.0
2-Methylnaphthalene	0.300	20.0	± 20.0	± 25.0
Acenaphthylene	0.900	20.0	± 20.0	± 25.0
Acenaphthene	0.500	20.0	± 20.0	± 25.0
Fluorene	0.700	20.0	± 25.0	± 50.0
Phenanthrene	0.300	20.0	± 25.0	± 50.0
Anthracene	0.400	20.0	± 25.0	± 50.0
Fluoranthene	0.400	20.0	± 25.0	± 50.0
Pyrene	0.500	20.0	± 30.0	± 50.0
Benzo(a)anthracene	0.400	20.0	± 25.0	± 50.0
Chrysene	0.400	20.0	± 25.0	± 50.0
Benzo(b)fluoranthene	0.100	20.0	± 30.0	± 50.0
Benzo(k)fluoranthene	0.100	20.0	± 30.0	± 50.0
Benzo(a)pyrene	0.100	20.0	± 25.0	± 50.0
Indeno(1,2,3-cd)pyrene	0.100	20.0	± 40.0	± 50.0
Dibenzo(a,h)anthracene	0.010	25.0	± 40.0	± 50.0
Benzo(g,h,i)perylene	0.020	25.0	± 40.0	± 50.0

DATA REVIEW WORKSHEETS

Pentachlorophenol	0.010	40.0	± 50.0	± 50.0
Deuterated Monitoring Compounds				
Analyte	Minimum RRF	Maximum %RSD	Opening Maximum %D ¹	Closing Maximum %D
1,4-Dioxane-d ₈	0.010	20.0	± 25.0	± 50.0
Phenol-d ₅	0.010	20.0	± 25.0	± 25.0
Bis-(2-chloroethyl)ether-d ₈	0.100	20.0	± 20.0	± 25.0
2-Chlorophenol-d ₄	0.200	20.0	± 20.0	± 25.0
4-Methylphenol-d ₈	0.010	20.0	± 20.0	± 25.0
4-Chloroaniline-d ₄	0.010	40.0	± 40.0	± 50.0
Nitrobenzene-d ₅	0.050	20.0	± 20.0	± 25.0
2-Nitrophenol-d ₄	0.050	20.0	± 20.0	± 25.0
2,4-Dichlorophenol-d ₃	0.060	20.0	± 20.0	± 25.0
Dimethylphthalate-d ₈	0.300	20.0	± 20.0	± 25.0
Acenaphthylene-d ₈	0.400	20.0	± 20.0	± 25.0
4-Nitrophenol-d ₄	0.010	40.0	± 40.0	± 50.0
Fluorene-d ₁₀	0.100	20.0	± 20.0	± 25.0
4,6-Dinitro-2-methylphenol-d ₂	0.010	40.0	± 30.0	± 50.0
Anthracene-d ₁₀	0.300	20.0	± 20.0	± 25.0
Pyrene-d ₁₀	0.300	20.0	± 25.0	± 50.0
Benzo(a)pyrene-d ₁₂	0.010	20.0	± 20.0	± 50.0
Fluoranthene-d ₁₀ (SIM)	0.400	20.0	± 25.0	± 50.0
2-Methylnaphthalene-d ₁₀ (SIM)	0.300	20.0	± 20.0	± 25.0

¹ If a closing CCV is acting as an opening CCV, all target analytes must meet the requirements for an opening CCV.

Note: If analysis by SIM technique is requested for PAH/pentachlorophenols, calibration standards analyzed at 0.10, 0.20, 0.40, 0.80, and 1.0 ng/uL for each target compound of interest and the associated DMCs. Pentachlorophenol will require only a four point initial calibration at 0.20, 0.40, 0.80, and 1.0 ng/uL.

DATA REVIEW WORKSHEETS

All criteria were met _____
 Criteria were not met _____
 and/or see below X

CONTINUING CALIBRATION VERIFICATION

Compliance requirements for satisfactory instrument calibration are established to ensure that the instrument is capable of producing and maintaining acceptable quantitative data.

Date of initial calibration: 04/19/16; 04/25/16 (Scan) 04/11/16 (Scan)
 Date of initial calibration verification (CCV): 04/19-20/16; 04/25/16 04/11/16
 Date of continuing calibration verification (CCV): 04/29/16 05/03/16
 Date of closing CCV: - -
 Instrument ID numbers: GCMS2P GCMS5P
 Matrix/Level: Aqueous/low Aqueous/low

Date of initial calibration: 04/12/16; 04/26/16 (Scan)
 Date of initial calibration verification (CCV): 04/12-13/16; 04/26-27/16
 Date of continuing calibration verification (CCV): 04/29/16; 05/06/16
 Date of closing CCV: -
 Instrument ID numbers: GCMSM
 Matrix/Level: Aqueous/low

Date of initial calibration: 04/21/16 (SIM) 04/14/16
 Date of initial calibration verification (CCV): 04/21/16 04/14/16
 Date of continuing calibration verification (CCV): 04/28/16 05/09/16; 05/10/16
 Date of closing CCV: -
 Instrument ID numbers: GCMS3M GCMS4M
 Matrix/Level: Aqueous/low Aqueous/low

DATE	LAB ID#	FILE	CRITERIA OUT RFs, %RSD, %D, r	COMPOUND	SAMPLES AFFECTED
See enclosed list					

Note: Initial and continuing calibration verifications meet the required criteria except the cases describe in the list enclosed. Analytes not detected in affected samples, results qualified (UJ).

No closing calibration verification included in data package. No action taken, professional judgment.

* Analytes with % difference in the continue calibration verification outside the method performance criteria but within the validation guidelines criteria, + 40 %. No action taken

DATA REVIEW WORKSHEETS

Actions:

Notes: Verify that the CCV is run at the required frequency (an opening and closing CCV must be run within 12-hour period).

All DMCs must meet the RRF values given in Table 2. No qualification of the data is necessary on DMCs RRF and %RSD/%D alone. Use professional judgment to evaluate DMCs and %RSD/%D data in conjunction with DMCs recoveries to determine the need for qualification of the data.

Qualify the initial calibration analytes listed in Table 2 using the following criteria in the CCVs:

Table 4. CCV Actions for Semivolatile Analysis

Criteria for Opening CCV	Criteria for Closing CCV	Action	
		Detect	Non-detect
CCV not performed at required frequency and sequence	CCV not performed at required frequency	Use professional judgment R	Use professional judgment R
CCV not performed at specified concentration	CCV not performed at specified concentration	Use professional judgment	Use professional judgment
RRF < Minimum RRF in Table 2 for target analyte	RRF < Minimum RRF in Table 2 for target analyte	Use professional judgment J or R	R
RRF ≥ Minimum RRF in Table 2 for target analyte	RRF ≥ Minimum RRF in Table 2 for target analyte	No qualification	No qualification
%D outside the Opening Maximum %D limits in Table 2 for target analyte	%D outside the Closing Maximum %D limits in Table 2 for target analyte	J	UJ
%D within the inclusive Opening Maximum %D limits in Table 2 for target analyte	%D within the inclusive Closing Maximum %D limits in Table 2 for target analyte	No qualification	No qualification

DATA REVIEW WORKSHEET

CONTINUING CALIBRATION VERIFICATION

INSTRUMENT: GCMS2P
DATE: 04/29/16
FILE ID: CC2547-50

Compound	%Dev
1,4-Dioxane*	24.2
Pyridine	23.2#
4-Nitrophenol*	-28.4#
4-Nitroaniline*	23.2#

CONTINUING CALIBRATION VERIFICATION

INSTRUMENT: GCMS5P
DATE: 05/03/16
FILE ID: cc1382-25

Compound	%Dev
N-Nitroso-di-n-propylamine	-22.5#
Caprolactam	-47.4#
2-Nitroaniline	-32.9#

CONTINUING CALIBRATION VERIFICATION

INSTRUMENT: GCMSM
DATE: 04/29/16
FILE ID: cc5217-25

Compound	%Dev
Hexachlorocyclopentadiene	-46.1#

INSTRUMENT: GCMSM
DATE: 04/29/16
FILE ID: cc5239-25

Compound	%Dev
----------	------

Atrazine*	-20.1#
-----------	--------

INSTRUMENT: GCMSM
DATE: 05/06/16
FILE ID: cc5239-50

Compound	%Dev
----------	------

Atrazine*	39.6#
-----------	-------

DATE: 05/06/16
FILE ID: cc5217-50

Compound	%Dev
----------	------

2,4-dinitrophenol*	-22.8
Pentachlorophenol*	-25.7#

DATA REVIEW WORKSHEETS

All criteria were met _____
 Criteria were not met _____
 and/or see below ____X____

BLANK ANALYSIS RESULTS (Sections 1 & 2)

The assessment of the blank analysis results is to determine the existence and magnitude of contamination problems. The criteria for evaluation of blanks apply only to blanks associated with the samples, including trip, equipment, and laboratory blanks. If problems with any blanks exist, all data associated with the case must be carefully evaluated to determine whether or not there is an inherent variability in the data for the case, or if the problem is an isolated occurrence not affecting other data.

List the contamination in the blanks below. High and low levels blanks must be treated separately.

Notes: The concentration of non-target compounds in all blanks must be less than or equal to 10 ug/L.

The concentration of target compounds in all blanks must be less than its CRQL listed in the method.

Samples taken from a drinking water tap do not have and associated field blank.

Laboratory blanks

DATE ANALYZED	LAB ID	LEVEL/ MATRIX	COMPOUND	CONCENTRATION UNITS
_No_target_analytes_detected_in_method_blanks_except_for_the_following:_____				
_04/29/16_____	_OP93453-MB1_____	_Aq/low_____	_Di-n-butyl_phthalate_____	_0.79_ug/L_____

Note: No action taken, concentration in blank < action level/reporting limits.

Field/Equipment/Trip blank

DATE ANALYZED	LAB ID	LEVEL/ MATRIX	COMPOUND	CONCENTRATION UNITS
_No_field/trip_blanks_analyzed_with_this_data_package._No_target_analyte_detected_in_the____ _equipment_blank_except_for_the_following:_____				
_04/29/16_____	_JC19023-3_____	_Aq./low_____	_Di-n-butyl_phthalate_____	_1.8_ug/L_____

Note: No action taken, concentration in blank < action level/reporting limits.

DATA REVIEW WORKSHEETS

All criteria were met _____
 Criteria were not met _____
 and/or see below X

BLANK ANALYSIS RESULTS (Section 3)

Blank Actions

Qualify samples based on the criteria summarized in Table 5:

Table 5. Blank and TCLP/SPLP LEB Actions for Semivolatile Analysis

Blank Type	Blank Result	Sample Result	Action
Method, TCLP/SPLP LEB, Field	Detect	Non-detect	No qualification
	< CRQL	< CRQL	Report at CRQL and qualify as non-detect (U)
		≥ CRQL	Use professional judgment
	≥ CRQL	< CRQL	Report at CRQL and qualify as non-detect (U)
		≥ CRQL but < Blank Result	Report at sample results and qualify as non-detect (U) or as unusable (R)
		≥ CRQL and ≥ Blank Result	Use professional judgment
	Grossly high	Detect	Report at sample results and qualify as unusable (R)
	TIC > 5.0 ug/L (water) or 0.0050 mg/L (TCLP leachate) or TIC > 170 ug/Kg (soil)	Detect	Use professional judgment

List samples qualified

CONTAMINATION SOURCE/LEVEL	COMPOUND	CONC/UNITS	AL/UNITS	SQL	AFFECTED SAMPLES

All criteria were met X
 Criteria were not met
 and/or see below

SURROGATE SPIKE RECOVERIES – DEUTERATED MONITORING COMPOUNDS (DMCs)

Laboratory performance of individual samples is established by evaluation of surrogate spike recoveries – deuterated monitoring compounds. All samples are spiked with surrogate compounds prior to sample analysis. The accuracy of the analysis is measured by the surrogate percent recovery. Since the effects of the sample matrix are frequently outside the control of the laboratory and may present relatively unique problems, the validation of data is frequently subjective and demands analytical experience and professional judgment.

Notes: Recoveries for DMCs in samples and blanks must be within the limits specified in Table 6.

The recovery limits for any of the compounds listed in Table 6 may be expanded at any time during the period of performance if USEPA determines that the limits are too restrictive.

If a DMC is not added in the samples and blanks or the concentrations of DMCs in the samples and blank not the specified, use professional judgment in qualifying the data.

Table 7. DMC Actions for Semivolatile Analysis

Criteria	Action	
	Detect	Non-detect
%R < 10% (excluding DMCs with 10% as a lower acceptance limit)	J-	R
10% ≤ %R (excluding DMCs with 10% as a lower acceptance limit) < Lower Acceptance Limit	J-	UJ
Lower Acceptance limit ≤ %R ≤ Upper Acceptance Limit	No qualification	No qualification
%R > Upper Acceptance Limit	J+	No qualification

List the percent recoveries (%Rs) which do not meet the criteria for DMCs (surrogate) recovery.

Matrix: Groundwater

SAMPLE ID	SURROGATE COMPOUND	ACTION
<u>DMCs meet the required criteria. Non-deuterated surrogates added to the samples</u>		
<u>within laboratory recovery limits.</u>		

Table 8. Semivolatile DMCs and the Associated Target Analytes

1,4-Dioxane-d₈ (DMC-1)	Phenol-d₅ (DMC-2)	Bis(2-Chloroethyl) ether-d₃ (DMC-3)
1,4-Dioxane	Benzaldehyde Phenol	Bis(2-chloroethyl)ether 2,2'-Oxybis(1-chloropropane) Bis(2-chloroethoxy)methane
2-Chlorophenol-d₄ (DMC-4)	4-Methylphenol-d₅ (DMC-5)	4-Chloroaniline-d₄ (DMC-6)
2-Chlorophenol	2-Methylphenol 3-Methylphenol 4-Methylphenol 2,4-Dimethylphenol	4-Chloroaniline Hexachlorocyclopentadiene Dichlorobenzidine
Nitrobenzene-d₅ (DMC-7)	2-Nitrophenol-d₄ (DMC-8)	2,4-Dichlorophenol-d₃ (DMC-9)
Acetophenone N-Nitroso-di-n-propylamine Hexachloroethane Nitrobenzene 2,6-Dinitrotoluene 2,4-Dinitrotoluene N-Nitrosodiphenylamine	Isophorone 2-Nitrophenol	2,4-Dichlorophenol Hexachlorobutadiene Hexachlorocyclopentadiene 4-Chloro-3-methylphenol 2,4,6-Trichlorophenol 2,4,5-Trichlorophenol 1,2,4,5-Tetrachlorobenzene *Pentachlorophenol 2,3,4,6-Tetrachlorophenol
Dimethylphthalate-d₆ (DMC-10)	Acenaphthylene-d₈ (DMC-11)	4-Nitrophenol-d₄ (DMC-12)
Caprolactam 1,1'-Biphenyl Dimethylphthalate Diethylphthalate Di-n-butylphthalate Butylbenzylphthalate Bis(2-ethylhexyl) phthalate Di-n-octylphthalate	*Naphthalene *2-Methylnaphthalene 2-Chloronaphthalene *Acenaphthylene *Acenaphthene	2-Nitroaniline 3-Nitroaniline 2,4-Dinitrophenol 4-Nitrophenol 4-Nitroaniline

DATA REVIEW WORKSHEETS

Fluorene-d₁₀ (DMC-13)	4,6-Dinitro-2-methylphenol-d₂ (DMC-14)	Anthracene-d₁₀ (DMC-15)
Dibenzofuran *Fluorene 4-Chlorophenyl-phenylether 4-Bromophenyl-phenylether Carbazole	4,6-Dinitro-2-methylphenol	Hexachlorobenzene Atrazine *Phenanthrene *Anthracene
Pyrene-d₁₀ (DMC-16)	Benzo(a)pyrene-d₁₂ (DMC-17)	
*Fluoranthene *Pyrene *Benzo(a)anthracene *Chrysene	3,3'-Dichlorobenzidine *Benzo(b)fluoranthene *Benzo(k)fluoranthene *Benzo(a)pyrene *Indeno(1,2,3-cd)pyrene *Dibenzo(a,h)anthracene *Benzo(g,h,i)perylene	

*Included in optional Target Analyte List (TAL) of PAHs and PCP only.

Table 9. Semivolatile SIM DMCs and the Associated Target Analytes

Fluoranthene-d₁₀ (DMC-1)	2-Methylnaphthalene-d₁₀ (DMC-2)
Fluoranthene	Naphthalene
Pyrene	2-Methylnaphthalene
Benzo(a)anthracene	Acenaphthylene
Chrysene	Acenaphthene
Benzo(b)fluoranthene	Fluorene
Benzo(k)fluoranthene	Pentachlorophenol
Benzo(a)pyrene	Phenanthrene
Indeno(1,2,3-cd)pyrene	Anthracene
Dibenzo(a,h)anthracene	
Benzo(g,h,i)perylene	

All criteria were met X
 Criteria were not met
 and/or see below _____

VII. A MATRIX SPIKE/MATRIX SPIKE DUPLICATE (MS/MSD)

This data is generated to determine long term precision and accuracy in the analytical method for various matrices. This data alone cannot be used to evaluate the precision and accuracy of individual samples. If any % R in the MS or MSD falls outside the designated range, the reviewer should determine if there are matrix effects, i.e. LCS data are within the QC limits but MS/MSD data are outside QC limit.

1. MS/MSD Recoveries and Precision Criteria

The laboratory should use one MS and a duplicate analysis of an unspiked field sample if target analytes are expected in the sample. If target analytes are not expected, MS/MSD should be analyzed.

NOTES: Data for MS and MSDs will not be present unless requested by the Region.
 Notify the Contract Laboratory COR if a field or trip blank was used for the MS and MSD.

For a Matrix Spike that does not meet criteria, apply the action to only the field sample used to prepare the Matrix Spike sample. If it is clearly stated in the data validation materials that the samples were taken through incremental sampling or some other method guaranteeing the homogeneity of the sample group, then the entire sample group may be qualified.

List the %Rs, RPD of the compounds which do not meet the criteria.

Sample ID: <u>JC19098-2</u>	Matrix/Level: <u>Aqueous</u>
Sample ID: <u>JC19138-2 (SIM)</u>	Matrix/Level: <u>Aqueous</u>
Sample ID: <u>JC18964-19</u>	Matrix/Level: <u>Soil</u>
Sample ID: <u>JC19023-1 (SIM)</u>	Matrix/Level: <u>Soil</u>

MS OR MSD	COMPOUND	% R	RPD	QC LIMITS	ACTION
<u>See enclosed list</u>					

Note: Results for 1,4-Dioxane qualified estimated (J) in sample JC19023-1.

- * QC limits are laboratory in-house performance criteria, LL = lower limit, UL = upper limit.
- * If QC limits are not available, use limits of 70 – 130 %.

Actions:

QUALITY	%R < LL	%R > UL
Positive results	J	J
Nondetects results	R	Accept

DATA REVIEW WORKSHEET

MATRIX SPIKE/MATRIX SPIKE DUPLICATE

SDG: JC19023

Matrix: Aqueous

Compound	JC19098-2 ug/l	Q	Spike ug/l	MS ug/l	MS %	Spike ug/l	MSD ug/l	MSD %	RPD	Limits Rec/RPD
Benzaldehyde	ND		100	330	330* a	100	452	452* a	31* b	10-164/30
Benzo(a)anthracene	ND		100	60.8	61	100	84.9	85	33* b	40-123/24
Benzo(a)pyrene	ND		100	61.9	62	100	83.5	84	30* b	41-127/25
Benzo(g,h,i)perylene	ND		100	51.4	51	100	76.1	76	39* b	34-128/28
Benzo(k)fluoranthene	ND		100	67.5	68	100	92.6	93	31* b	39-122/26
Butyl benzyl phthalate	ND		100	69.3	69	100	94.6	95	31* b	21-146/28
Chrysene	ND		100	59.6	60	100	82.9	83	33* b	41-128/24
Dibenzo(a,h)- anthracene	ND		100	57.4	57	100	81.8	82	35* b	35-130/27
bis(2-Ethylhexyl)- phthalate	ND		100	69.8	70	100	95.5	96	31* b	34-141/28
Indeno(1,2,3-cd)pyrene	ND		100	54.1	54	100	79.3	79	38* b	32-130/30
1-Methylnaphthalene	ND		100	139	139* a	100	165	165* a	17	34-124/25
2-Methylnaphthalene	128 c		100	201	72	100	253	124* a	23	34-123/24
Pyrene	ND		100	64.7	65	100	89.7	90	32* b	43-124/26

* = Outside of Control Limits.

(a) Outside control limits due to matrix interference.

(b) Outside of in house control limits.

(c) Result is from Run #2.

SDG: JC19023

Matrix: Soil

Compound	JC18964-19 ug/kg	Q	Spike ug/kg	MS ug/kg	MS %	Spike ug/kg	MSD ug/kg	MSD %	RPD	Limits Rec/RPD
3&4-Methylphenol	4350 b		1650	6850	151* a	1690	8100	222* a	17	32-113/34
Phenol	1300		1650	4270	180* c	1690	5070	224* c	17	25-112/33
bis(2-Chloroethyl)ether	ND		1650	1550	94	1690	2130	126* c	32	30-113/37
3,3'-Dichlorobenzidine	ND		3310	239	7* c	3370	303	9* c	24	10-115/44

(a) Outside control limits due to high level in sample relative to spike amount.

(b) Result is from Run #2.

(c) Outside control limits due to matrix interference.

The QC reported here applies to the following samples:
JC19023-3

Method: **SW846 8270D BY SIM**

SDG: JC19023

Matrix: Aqueous

Compound	JC19138-3 ug/l	Q	Spike ug/l	MS ug/l	MS %	Spike ug/l	MSD ug/l	MSD %	RPD	Limits Rec/RPD
Naphthalene	1.35		2	0.926	0* a	2	1.27	0* a	31	23-140/36
1,4-Dioxane	ND		2	0.609	30	2	0.992	50	48* b	20-160/30

The QC reported here applies to the following samples:
JC19023-1, JC19023-2, JC19023-4

Method: **SW846 8270D BY SIM**

Compound	JC19023-1 ug/kg	Q	Spike ug/kg	MS ug/kg	MS %	Spike ug/kg	MSD ug/kg	MSD %	RPD	Limits Rec/RPD
1,4-Dioxane	7.80		213	19.9	6* a	213	25.8	8* a	26	50-150/30
Naphthalene	30.0		42.6	40.9	26	42.6	40.8	25	0	10-190/36

- (a) Outside control limits due to matrix interference.
(b) Analytical precision exceeds in-house control limits.

DATA REVIEW WORKSHEETS

MS/MSD criteria apply only to the unspiked sample, its dilutions, and the associated MS/MSD samples:

If the % R for the affected compounds were < LL (or 70 %), qualify positive results (J) and nondetects (UJ).

If the % R for the affected compounds were > UL (or 130 %), only qualify positive results (J).

If 25 % or more of all MS/MSD %R were < LL (or 70 %) or if two or more MS/MSD %Rs were < 10%, qualify all positive results (J) and reject nondetects (R).

A separate worksheet should be used for each MS/MSD pair.

All criteria were met ☒ X
Criteria were not met
and/or see below _____

INTERNAL STANDARD PERFORMANCE

The assessment of the internal standard (IS) parameter is used to assist the data reviewer in determining the condition of the analytical instrumentation.

List the internal standard area of samples which do not meet the criteria.

DATE	SAMPLE ID	IS OUT	IS AREA	ACCEPTABLE RANGE	ACTION
------	-----------	--------	---------	------------------	--------

Internal standard area counts meet the required criteria.

Action:

1. If an internal standard area count for a sample or blank is greater than 200.0% of the area for the associated standard (opening CCV or mid-point standard from initial calibration) (see Table 10 below):
 - a. Qualify detects for compounds quantitated using that internal standard as estimated low (J-).
 - b. Do not qualify non-detected associated compounds.
2. If an internal standard area count for a sample or blank is less than 20.0% of the area for the associated standard (opening CCV or mid-point standard from initial calibration):
 - a. Qualify detects for compounds quantitated using that internal standard as estimated high (J+).
 - b. Qualify non-detected associated compounds as unusable (R).
3. If an internal standard area count for a sample or blank is greater than or equal to 50.0%, and less than or equal to 200% of the area for the associated standard opening CCV or mid-point standard from initial calibration, no qualification of the data is necessary.
4. If an internal standard RT varies by more than 10.0 seconds: Examine the chromatographic profile for that sample to determine if any false positives or negatives exist. For shifts of a large magnitude, the reviewer may consider partial or total rejection of the data for that sample fraction. Detects should not need to be qualified as unusable (R) if the mass spectral criteria are met.
5. If an internal standard RT varies by less than or equal to 10.0 seconds, no qualification of the data is necessary.

Note: Inform the Contract Laboratory Program Project Officer (CLP PO) if the internal standard performance criteria are grossly exceeded. Note in the Data Review Narrative potential effects on the data resulting from unacceptable internal standard performance.

DATA REVIEW WORKSHEETS

State in the Data Review Narrative if the required internal standard compounds are not added to a sample or blank or if the required internal standard compound is not analyzed at the specified concentration.

Actions:

Table 10. Internal Standard Actions for Semivolatile Analysis

Criteria	Action	
	Detect	Non-detect
Area response < 20% of the opening CCV or mid-point standard CS3 from ICAL	J+	R
20% ≤ Area response < 50% of the opening CCV or mid-point standard CS3 from ICAL	J+	UJ
50% ≤ Area response ≤ 200% of the opening CCV or mid-point standard CS3 from ICAL	No qualification	No qualification
Area response > 200% of the opening CCV or mid-point standard CS3 from ICAL	J-	No qualification
RT shift between sample/blank and opening CCV or mid-point standard CS3 from ICAL > 10.0 seconds	R	R
RT shift between sample/blank and opening CCV or mid-point standard CS3 from ICAL < 10.0 seconds	No qualification	No qualification

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below _____

TARGET COMPOUND IDENTIFICATION

Criteria:

Is the Relative Retention Times (RRTs) of reported compounds within ± 0.06 RRT units of the standard RRT [opening Continuing Calibration Verification (CCV) or mid-point standard from the initial calibration].
Yes? or No?

List compounds not meeting the criteria described above:

Sample ID	Compounds	Actions
=====	=====	=====
_____	_____	_____
_____	_____	_____
_____	_____	_____

Mass spectra of the sample compound and a current laboratory-generated standard [i.e., the mass spectrum from the associated calibration standard (opening CCV or mid-point standard from initial calibration)] must match according to the following criteria:

- All ions present in the standard mass spectrum at a relative intensity greater than 10% must be present in the sample spectrum.
- The relative intensities of these ions must agree within $\pm 20\%$ between the standard and sample spectra (e.g., for an ion with an abundance of 50% in the standard spectrum, the corresponding sample ion abundance must be between 30-70%).
- Ions present at greater than 10% in the sample mass spectrum, but not present in the standard spectrum, must be evaluated by a reviewer experienced in mass spectral interpretation.

List compounds not meeting the criteria described above:

Sample ID	Compounds	Actions
=====	=====	=====
_____	_____	_____
_____	_____	_____
_____	_____	_____

Identified compounds meet the required criteria

DATA REVIEW WORKSHEETS

Action:

1. The application of qualitative criteria for GC/MS analysis of target compounds requires professional judgment. It is up to the reviewer's discretion to obtain additional information from the laboratory. If it is determined that incorrect identifications were made, qualify all such data as unusable (R).
2. Use professional judgment to qualify the data if it is determined that cross-contamination has occurred.
3. Note in the Data Review Narrative any changes made to the reported compounds or concerns regarding target compound identifications. Note, for Contract Laboratory COR action, the necessity for numerous or significant changes.

TENTATIVELY IDENTIFIED COMPOUNDS (TICS)

NOTE: Tentatively identified compounds should only be evaluated when requested by a party from outside of the Hazardous Waste Support Section (HWSS).

List TICs

Sample ID	Compound	Sample ID	Compound
=====		=====	
_____		_____	
_____		_____	
_____		_____	
_____		_____	

Action:

1. Qualify all TIC results for which there is presumptive evidence of a match (e.g. greater than or equal to 85% match) as tentatively identified (NJ), with approximated concentrations. TICs labeled "unknown" are qualified as estimated (J).
2. General actions related to the review of TIC results are as follows:
 - a. If it is determined that a tentative identification of a non-target compound is unacceptable, change the tentative identification to "unknown" or another appropriate identification, and qualify the result as estimated (J).
 - b. If all contractually-required peaks were not library searched and quantitated, the Region's designated representative may request these data from the laboratory.
3. In deciding whether a library search result for a TIC represents a reasonable identification, use professional judgment. If there is more than one possible match, report the result as "either compound X or compound Y". If there is a lack of isomer specificity, change the TIC result to a nonspecific isomer result (e.g., 1,3,5-trimethyl benzene to trimethyl benzene isomer) or to a compound class (e.g., 2-methyl, 3-ethyl benzene to a substituted aromatic compound).
4. The reviewer may elect to report all similar compounds as a total (e.g., all alkanes may be summarized and reported as total hydrocarbons).

DATA REVIEW WORKSHEETS

5. Target compounds from other fractions and suspected laboratory contaminants should be marked as "non-reportable".
6. Other Case factors may influence TIC judgments. If a sample TIC match is poor, but other samples have a TIC with a valid library match, similar RRT, and the same ions, infer identification information from the other sample TIC results.
7. Note in the Data Review Narrative any changes made to the reported data or any concerns regarding TIC identifications.
8. Note, for Contract Laboratory COR action, failure to properly evaluate and report TICs

All criteria were met X
 Criteria were not met
 and/or see below

SAMPLE QUANTITATION AND REPORTED CONTRACT REQUIRED QUANTITATION LIMITS (CRQLS)

Action:

1. When a sample is analyzed at more than one dilution, the lower CRQL are used unless a QC exceedance dictates the use of higher CRQLs from the diluted sample. Samples reported with an "E" qualifier should be reported from the diluted sample.
2. If any discrepancies are found, the Region's designated representative may contact the laboratory to obtain additional information that could resolve any differences. If a discrepancy remains unresolved, the reviewer must use professional judgment to decide which value is the most accurate. Under these circumstances, the reviewer may determine that qualification of data is warranted. Note in the Data Review Narrative a description of the reasons for data qualification and the qualification that is applied to the data.
3. For non-aqueous samples, if the solids is less than 10.0%, use professional judgment for both detects and non-detects. If the percent solid for a soil sample is greater than or equal to 10.0% and less than 30.0%, use professional judgment to qualify detects and non-detects. If the percent solid for a soil sample is greater than or equal to 30.0%, detects and non-detects should not be qualified (see Table 11).
4. Note, for Contract Laboratory COR action, numerous or significant failures to accurately quantify the target compounds or to properly evaluate and adjust CRQLs.
5. Results between MDL and CRQL should be qualified as estimated "J".
6. Results < MDL should be reported at the CRQL and qualified "U". MDLs themselves should not be reported.

Table 11. Percent Solids Actions for Semivolatile Analysis for Non-Aqueous Samples

Criteria	Action	
	Detects	Non-detects
%Solids < 10.0%	Use professional judgment	Use professional judgment
10.0% ≤ %Solids ≤ 30.0%	Use professional judgment	Use professional judgment
%Solids > 30.0%	No qualification	No qualification

SAMPLE QUANTITATION

The sample quantitation evaluation is to verify laboratory quantitation results. In the space below, please show a minimum of one sample calculation:

Sample ID: JC19023-1 Analyte: Benzaldehyde RF: 0.904

$$\begin{aligned}
 [] &= (21467)(40)/(404573)(0.904) \\
 &= 2.35 \text{ ppm} \quad \text{Ok}
 \end{aligned}$$

DATA REVIEW WORKSHEETS

QUANTITATION LIMITS

A. Dilution performed

[illegible]

DATA REVIEW WORKSHEETS

All criteria were met N/A
 Criteria were not met
 and/or see below

FIELD DUPLICATE PRECISION

Sample IDs: -

Matrix: -

Field duplicate samples may be taken and analyzed as an indication of overall precision. These analyses measure both field and lab precision; therefore, the results may have more variability than laboratory duplicates which only laboratory performance. It is also expected that soil duplicate results will have a greater variance than water matrices due to difficulties associated with collecting identical field duplicate samples.

The project QAPP should be reviewed for project-specific information.

Suggested criteria: if large RPD (> 50 %) is observed, confirm identification of the samples and note differences. If both samples and duplicate are <5 SQL, the RPD criteria is doubled.

COMPOUND	SQL ug/L	SAMPLE CONC.	DUPLICATE CONC.	RPD	ACTION
No field/laboratory duplicate analyzed as part of this data package. MS/MSD % recoveries RPD used to assess precision; RPD within the required criteria < 50 % for detected target analytes.					

DATA REVIEW WORKSHEETS

All criteria were met ☒
 Criteria were not met
 and/or see below _____

OTHER ISSUES

A. System Performance

List samples qualified based on the degradation of system performance during sample analysis:

Sample ID	Comments	Actions
=====	=====	=====
_____	_____	_____
_____	_____	_____
_____	_____	_____
_____	_____	_____

Action:

Use professional judgment to qualify the data if it is determined that system performance has degraded during sample analyses. Inform the Contract Laboratory Program COR any action as a result of degradation of system performance which significantly affected the data.

B. Overall Assessment of Data

List samples qualified based on other issues:

Sample ID	Comments	Actions
=====	=====	=====
_____	_____	_____
_____	_____	_____
_____	_____	_____

Action:

1. Use professional judgment to determine if there is any need to qualify data which were not qualified based on the Quality Control (QC) criteria previously discussed.
2. Write a brief narrative to give the user an indication of the analytical limitations of the data. Inform the Contract Laboratory COR the action, any inconsistency of the data with the Sample Delivery Group (SDG) Narrative. If sufficient information on the intended use and required quality of the data is available, the reviewer should include their assessment of the usability of the data within the given context. This may be used as part of a formal Data Quality Assessment (DQA).
3. Sometimes, due to dilutions, re-analysis or SIM/Scan runs are being performed, there will be multiple results for a single analyte from a single sample. The following criteria and professional judgment are used to determine which result should be reported:
 - The analysis with the lower CRQL
 - The analysis with the better QC results
 - The analysis with the higher results