

CERTIFICATION

SDG No: JC19473 Laboratory: Accutest, New Jersey
 Site: BMS, Building 5 Area, PR Matrix: Accutest, Florida
 Humacao, PR Soil/Groundwater

SUMMARY: Groundwater and soil samples (Table 1) were collected on the BMSMC facility – Building 5 Area. The BMSMC facility is located in Humacao, PR. Samples were taken April 30-May 2, 2016 and were analyzed in Accutest Laboratory of Dayton, New Jersey for the ABN TCL Special List (1,4-Dioxane and Naphthalene were analyzed following the SIM technique) and for TCL pesticides list that reported the data under SDG No.: JC19473. Accutest, Florida analyzed for low molecular weight alcohols (LMWA), results were also reported under SDG No.: JC19473. Results were validated using the latest validation guidelines (July, 2015) of the EPA Hazardous Waste Support Section. The analyses performed are shown in Table 1. Individual data review worksheets are enclosed for each target analyte group. The data sample organic data samples summary form shows for analytes results that were qualified.

In summary the results are valid and can be used for decision taking purposes.

Table 1. Samples analyzed and analysis performed

SAMPLE ID	SAMPLE DESCRIPTION	MATRIX	ANALYSIS PERFORMED
JC19473-1	BPEB-10	AQ – Equipment blank	ABN TCL special list; pesticides TCL list; LMWA
JC19473-2	RA-18 GWS	Groundwater	ABN TCL special list; pesticides TCL list;
JC19473-2A	RA-18 GWS	Groundwater	LMWA
JC19473-3	RA-18 (12.5-13.5)	Soil	ABN TCL special list; pesticides TCL list; LMWA

Reviewer Name: Rafael Infante
 Chemist License 1888

Signature:

Date:

Rafael Infante

May 23, 2016



SGS Accutest

Report of Analysis

Page 1 of 3

Client Sample ID: BPEB-10
 Lab Sample ID: JC19473-1
 Matrix: AQ - Equipment Blank
 Method: SW846 8270D SW846 3510C
 Project: BMSMC, Building 5 Area, PR

Date Sampled: 04/30/16
 Date Received: 05/03/16
 Percent Solids: n/a

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	F156757.D	1	05/03/16	LK	05/03/16	OP93555	EF6600
Run #2							

Run #	Initial Volume	Final Volume
Run #1	900 ml	1.0 ml
Run #2		

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units	Q
95-57-8	2-Chlorophenol	ND	5.6	0.91	ug/l	
59-50-7	4-Chloro-3-methyl phenol	ND	5.6	0.99	ug/l	
120-83-2	2,4-Dichlorophenol	ND	2.2	1.4	ug/l	
105-67-9	2,4-Dimethylphenol	ND	5.6	2.7	ug/l	
51-28-5	2,4-Dinitrophenol	ND	11	1.7	ug/l	
534-52-1	4,6-Dinitro-o-cresol	ND	5.6	1.4	ug/l	
95-48-7	2-Methylphenol	ND	2.2	0.99	ug/l	
	3&4-Methylphenol	ND	2.2	0.98	ug/l	
88-75-5	2-Nitrophenol	ND	5.6	1.1	ug/l	
100-02-7	4-Nitrophenol	ND	11	1.3	ug/l	
87-86-5	Pentachlorophenol	ND	5.6	1.5	ug/l	
108-95-2	Phenol	ND	2.2	0.44	ug/l	
58-90-2	2,3,4,6-Tetrachlorophenol	ND	5.6	1.6	ug/l	
95-95-4	2,4,5-Trichlorophenol	ND	5.6	1.5	ug/l	
88-06-2	2,4,6-Trichlorophenol	ND	5.6	1.0	ug/l	
83-32-9	Acenaphthene	ND	1.1	0.21	ug/l	
208-96-8	Acenaphthylene	ND	1.1	0.15	ug/l	
98-86-2	Acetophenone	ND	2.2	0.23	ug/l	
120-12-7	Anthracene	ND	1.1	0.23	ug/l	
1912-24-9	Atrazine	ND	2.2	0.50	ug/l	
100-52-7	Benzaldehyde	ND	5.6	0.32	ug/l	
56-55-3	Benzo(a)anthracene	ND	1.1	0.23	ug/l	
50-32-8	Benzo(a)pyrene	ND	1.1	0.24	ug/l	
205-99-2	Benzo(b)fluoranthene	ND	1.1	0.23	ug/l	
191-24-2	Benzo(g,h,i)perylene	ND	1.1	0.38	ug/l	
207-08-9	Benzo(k)fluoranthene	ND	1.1	0.23	ug/l	
101-55-3	4-Bromophenyl phenyl ether	ND	2.2	0.45	ug/l	
85-68-7	Butyl benzyl phthalate	ND	2.2	0.51	ug/l	
92-52-4	1,1'-Biphenyl	ND	1.1	0.24	ug/l	
91-58-7	2-Chloronaphthalene	ND	2.2	0.26	ug/l	
106-47-8	4-Chloroaniline	ND	5.6	0.38	ug/l	
86-74-8	Carbazole	ND	1.1	0.25	ug/l	



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

Report of Analysis

Client Sample ID: BPEB-10
 Lab Sample ID: JC19473-1
 Matrix: AQ - Equipment Blank
 Method: SW846 8270D SW846 3510C
 Project: BSMC, Building 5 Area, PR

Date Sampled: 04/30/16
 Date Received: 05/03/16
 Percent Solids: n/a

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units	Q
105-60-2	Caprolactam	ND	2.2	0.72	ug/l	
218-01-9	Chrysene	ND	1.1	0.20	ug/l	
111-91-1	bis(2-Chloroethoxy)methane	ND	2.2	0.31	ug/l	
111-44-4	bis(2-Chloroethyl)ether	ND	2.2	0.28	ug/l	
108-60-1	bis(2-Chloroisopropyl)ether	ND	2.2	0.45	ug/l	
7005-72-3	4-Chlorophenyl phenyl ether	ND	2.2	0.41	ug/l	
121-14-2	2,4-Dinitrotoluene	ND	1.1	0.61	ug/l	
606-20-2	2,6-Dinitrotoluene	ND	1.1	0.53	ug/l	
91-94-1	3,3'-Dichlorobenzidine	ND	2.2	0.56	ug/l	
53-70-3	Dibenzo(a,h)anthracene	ND	1.1	0.37	ug/l	
132-64-9	Dibenzofuran	ND	5.6	0.24	ug/l	
84-74-2	Di-n-butyl phthalate	ND	2.2	0.55	ug/l	
117-84-0	Di-n-octyl phthalate	ND	2.2	0.26	ug/l	
84-66-2	Diethyl phthalate	ND	2.2	0.29	ug/l	
131-11-3	Dimethyl phthalate	ND	2.2	0.24	ug/l	
117-81-7	bis(2-Ethylhexyl)phthalate	ND	2.2	1.8	ug/l	
206-44-0	Fluoranthene	ND	1.1	0.19	ug/l	
86-73-7	Fluorene	ND	1.1	0.19	ug/l	
118-74-1	Hexachlorobenzene	ND	1.1	0.36	ug/l	
87-68-3	Hexachlorobutadiene	ND	1.1	0.55	ug/l	
77-47-4	Hexachlorocyclopentadiene	ND	11	3.1	ug/l	
67-72-1	Hexachloroethane	ND	2.2	0.43	ug/l	
193-39-5	Indeno(1,2,3-cd)pyrene	ND	1.1	0.37	ug/l	
78-59-1	Isophorone	ND	2.2	0.31	ug/l	
90-12-0	1-Methylnaphthalene	ND	1.1	0.29	ug/l	
91-57-6	2-Methylnaphthalene	ND	1.1	0.23	ug/l	
88-74-4	2-Nitroaniline	ND	5.6	0.31	ug/l	
99-09-2	3-Nitroaniline	ND	5.6	0.43	ug/l	
100-01-6	4-Nitroaniline	ND	5.6	0.49	ug/l	
98-95-3	Nitrobenzene	ND	2.2	0.71	ug/l	
621-64-7	N-Nitroso-di-n-propylamine	ND	2.2	0.53	ug/l	
86-30-6	N-Nitrosodiphenylamine	ND	5.6	0.25	ug/l	
85-01-8	Phenanthrene	ND	1.1	0.19	ug/l	
129-00-0	Pyrene	ND	1.1	0.24	ug/l	
95-94-3	1,2,4,5-Tetrachlorobenzene	ND	2.2	0.41	ug/l	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
367-12-4	2-Fluorophenol	65%		14-88%
4165-62-2	Phenol-d5	42%		10-110%

ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound



Report of Analysis

Client Sample ID:	BPEB-10	Date Sampled:	04/30/16
Lab Sample ID:	JC19473-1	Date Received:	05/03/16
Matrix:	AQ - Equipment Blank	Percent Solids:	n/a
Method:	SW846 8270D SW846 3510C		
Project:	BMSMC, Building 5 Area, PR		

ABN TCL Special List

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
118-79-6	2,4,6-Tribromophenol	114%		39-149%
4165-60-0	Nitrobenzene-d5	115%		32-128%
321-60-8	2-Fluorobiphenyl	97%		35-119%
1718-51-0	Terphenyl-d14	104%		10-126%



ND = Not detected MDL = Method Detection Limit
RL = Reporting Limit
E = Indicates value exceeds calibration range

J = Indicates an estimated value
B = Indicates analyte found in associated method blank
N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID:	BPEB-10	Date Sampled:	04/30/16
Lab Sample ID:	JC19473-1	Date Received:	05/03/16
Matrix:	AQ - Equipment Blank	Percent Solids:	n/a
Method:	SW846 8270D BY SIM SW846 3510C		
Project:	BMSMC, Building 5 Area, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	3M61049.D	1	05/04/16	LK	05/03/16	OP93555A	E3M2870
Run #2							

Run #	Initial Volume	Final Volume
Run #1	900 ml	1.0 ml
Run #2		

CAS No.	Compound	Result	RL	MDL	Units	Q
91-20-3	Naphthalene	ND	0.11	0.033	ug/l	
123-91-1	1,4-Dioxane	ND	0.11	0.054	ug/l	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
4165-60-0	Nitrobenzene-d5	92%		24-125%
321-60-8	2-Fluorobiphenyl	105%		19-127%
1718-51-0	Terphenyl-d14	115%		10-119%



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID:	BPEB-10	Date Sampled:	04/30/16
Lab Sample ID:	JC19473-1	Date Received:	05/03/16
Matrix:	AQ - Equipment Blank	Percent Solids:	n/a
Method:	SW846 8015C		
Project:	BMSMC, Building 5 Area, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1 ^a	XY064114.D	1	05/05/16	AFL	n/a	n/a	F:GXY2774
Run #2							

CAS No.	Compound	Result	RL	MDL	Units	Q
64-17-5	Ethanol	ND	5.0	1.0	mg/l	
78-83-1	Isobutyl Alcohol	ND	5.0	1.0	mg/l	
67-63-0	Isopropyl Alcohol	ND	5.0	1.0	mg/l	
71-23-8	n-Propyl Alcohol	ND	5.0	1.0	mg/l	
71-36-3	n-Butyl Alcohol	ND	5.0	1.0	mg/l	
67-56-1	Methanol	ND	5.0	1.0	mg/l	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
111-27-3	Hexanol	91%		73-123%

(a) Analysis performed at Accutest Laboratories, Orlando FL.



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID:	BPEB-10	Date Sampled:	04/30/16
Lab Sample ID:	JC19473-1	Date Received:	05/03/16
Matrix:	AQ - Equipment Blank	Percent Solids:	n/a
Method:	SW846 8081B SW846 3510C		
Project:	BMSMC, Building 5 Area, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	6G34643.D	1	05/03/16	DS	05/03/16	OP93571	G6G1001
Run #2							

Run #	Initial Volume	Final Volume
Run #1	300 ml	2.0 ml
Run #2		

Pesticide TCL List

CAS No.	Compound	Result	RL	MDL	Units	Q
309-00-2	Aldrin	ND	0.0067	0.0040	ug/l	
319-84-6	alpha-BHC	ND	0.0067	0.0040	ug/l	
319-85-7	beta-BHC	ND	0.0067	0.0038	ug/l	
319-86-8	delta-BHC	ND	0.0067	0.0030	ug/l	
58-89-9	gamma-BHC (Lindane)	ND	0.0067	0.0019	ug/l	
5103-71-9	alpha-Chlordane	ND	0.0067	0.0031	ug/l	
5103-74-2	gamma-Chlordane	ND	0.0067	0.0031	ug/l	
60-57-1	Dieldrin	ND	0.0067	0.0024	ug/l	
72-54-8	4,4'-DDD	ND	0.0067	0.0025	ug/l	
72-55-9	4,4'-DDE	ND	0.0067	0.0041	ug/l	
50-29-3	4,4'-DDT	ND	0.0067	0.0033	ug/l	
72-20-8	Endrin	ND	0.0067	0.0034	ug/l	
1031-07-8	Endosulfan sulfate	ND	0.0067	0.0035	ug/l	
7421-93-4	Endrin aldehyde	ND	0.0067	0.0034	ug/l	
53494-70-5	Endrin ketone	ND	0.0067	0.0034	ug/l	
959-98-8	Endosulfan-I	ND	0.0067	0.0033	ug/l	
33213-65-9	Endosulfan-II	ND	0.0067	0.0029	ug/l	
76-44-8	Heptachlor	ND	0.0067	0.0025	ug/l	
1024-57-3	Heptachlor epoxide	ND	0.0067	0.0044	ug/l	
72-43-5	Methoxychlor	ND	0.013	0.0038	ug/l	
8001-35-2	Toxaphene	ND	0.17	0.12	ug/l	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
877-09-8	Tetrachloro-m-xylene	57%		26-132%
877-09-8	Tetrachloro-m-xylene	58%		26-132%
2051-24-3	Decachlorobiphenyl	54%		10-118%
2051-24-3	Decachlorobiphenyl	50%		10-118%



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 3

Client Sample ID: RA-18 GWS
 Lab Sample ID: JC19473-2
 Matrix: AQ - Ground Water
 Method: SW846 8270D SW846 3510C
 Project: BSMC, Building 5 Area, PR

Date Sampled: 05/02/16
 Date Received: 05/03/16
 Percent Solids: n/a

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	F156758.D	1	05/03/16	LK	05/03/16	OP93555	EF6600
Run #2							

Run #	Initial Volume	Final Volume
Run #1	900 ml	1.0 ml
Run #2		

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units	Q
95-57-8	2-Chlorophenol	ND	5.6	0.91	ug/l	
59-50-7	4-Chloro-3-methyl phenol	ND	5.6	0.99	ug/l	
120-83-2	2,4-Dichlorophenol	ND	2.2	1.4	ug/l	
105-67-9	2,4-Dimethylphenol	ND	5.6	2.7	ug/l	
51-28-5	2,4-Dinitrophenol	ND	11	1.7	ug/l	
534-52-1	4,6-Dinitro-o-cresol	ND	5.6	1.4	ug/l	
95-48-7	2-Methylphenol	ND	2.2	0.99	ug/l	
	3&4-Methylphenol	ND	2.2	0.98	ug/l	
88-75-5	2-Nitrophenol	ND	5.6	1.1	ug/l	
100-02-7	4-Nitrophenol	ND	11	1.3	ug/l	
87-86-5	Pentachlorophenol	ND	5.6	1.5	ug/l	
108-95-2	Phenol	ND	2.2	0.44	ug/l	
58-90-2	2,3,4,6-Tetrachlorophenol	ND	5.6	1.6	ug/l	
95-95-4	2,4,5-Trichlorophenol	ND	5.6	1.5	ug/l	
88-06-2	2,4,6-Trichlorophenol	ND	5.6	1.0	ug/l	
83-32-9	Acenaphthene	ND	1.1	0.21	ug/l	
208-96-8	Acenaphthylene	ND	1.1	0.15	ug/l	
98-86-2	Acetophenone	ND	2.2	0.23	ug/l	
120-12-7	Anthracene	ND	1.1	0.23	ug/l	
1912-24-9	Atrazine	ND	2.2	0.50	ug/l	
100-52-7	Benzaldehyde	ND	5.6	0.32	ug/l	
56-55-3	Benzo(a)anthracene	ND	1.1	0.23	ug/l	
50-32-8	Benzo(a)pyrene	ND	1.1	0.24	ug/l	
205-99-2	Benzo(b)fluoranthene	ND	1.1	0.23	ug/l	
191-24-2	Benzo(g,h,i)perylene	ND	1.1	0.38	ug/l	
207-08-9	Benzo(k)fluoranthene	ND	1.1	0.23	ug/l	
101-55-3	4-Bromophenyl phenyl ether	ND	2.2	0.45	ug/l	
85-68-7	Butyl benzyl phthalate	ND	2.2	0.51	ug/l	
92-52-4	1,1'-Biphenyl	ND	1.1	0.24	ug/l	
91-58-7	2-Chloronaphthalene	ND	2.2	0.26	ug/l	
106-47-8	4-Chloroaniline	ND	5.6	0.38	ug/l	
86-74-8	Carbazole	ND	1.1	0.25	ug/l	

ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound



Report of Analysis

Client Sample ID: RA-18 GWS
 Lab Sample ID: JC19473-2
 Matrix: AQ - Ground Water
 Method: SW846 8270D SW846 3510C
 Project: BMSMC, Building 5 Area, PR

Date Sampled: 05/02/16
 Date Received: 05/03/16
 Percent Solids: n/a

AEN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units	Q
105-60-2	Caprolactam	ND	2.2	0.72	ug/l	
218-01-9	Chrysene	ND	1.1	0.20	ug/l	
111-91-1	bis(2-Chloroethoxy)methane	ND	2.2	0.31	ug/l	
111-44-4	bis(2-Chloroethyl)ether	ND	2.2	0.28	ug/l	
108-60-1	bis(2-Chloroisopropyl)ether	ND	2.2	0.45	ug/l	
7005-72-3	4-Chlorophenyl phenyl ether	ND	2.2	0.41	ug/l	
121-14-2	2,4-Dinitrotoluene	ND	1.1	0.61	ug/l	
606-20-2	2,6-Dinitrotoluene	ND	1.1	0.53	ug/l	
91-94-1	3,3'-Dichlorobenzidine	ND	2.2	0.56	ug/l	
53-70-3	Dibenzo(a,h)anthracene	ND	1.1	0.37	ug/l	
132-64-9	Dibenzofuran	ND	5.6	0.24	ug/l	
84-74-2	Di-n-butyl phthalate	ND	2.2	0.55	ug/l	
117-84-0	Di-n-octyl phthalate	ND	2.2	0.26	ug/l	
84-66-2	Diethyl phthalate	ND	2.2	0.29	ug/l	
131-11-3	Dimethyl phthalate	ND	2.2	0.24	ug/l	
117-81-7	bis(2-Ethylhexyl)phthalate	ND	2.2	1.8	ug/l	
206-44-0	Fluoranthene	ND	1.1	0.19	ug/l	
86-73-7	Fluorene	ND	1.1	0.19	ug/l	
118-74-1	Hexachlorobenzene	ND	1.1	0.36	ug/l	
87-68-3	Hexachlorobutadiene	ND	1.1	0.55	ug/l	
77-47-4	Hexachlorocyclopentadiene	ND	11	3.1	ug/l	
67-72-1	Hexachloroethane	ND	2.2	0.43	ug/l	
193-39-5	Indeno(1,2,3-cd)pyrene	ND	1.1	0.37	ug/l	
78-59-1	Isophorone	ND	2.2	0.31	ug/l	
90-12-0	1-Methylnaphthalene	ND	1.1	0.29	ug/l	
91-57-6	2-Methylnaphthalene	ND	1.1	0.23	ug/l	
88-74-4	2-Nitroaniline	ND	5.6	0.31	ug/l	
99-09-2	3-Nitroaniline	ND	5.6	0.43	ug/l	
100-01-6	4-Nitroaniline	ND	5.6	0.49	ug/l	
98-95-3	Nitrobenzene	ND	2.2	0.71	ug/l	
621-64-7	N-Nitroso-di-n-propylamine	ND	2.2	0.53	ug/l	
86-30-6	N-Nitrosodiphenylamine	ND	5.6	0.25	ug/l	
85-01-8	Phenanthrene	ND	1.1	0.19	ug/l	
129-00-0	Pyrene	ND	1.1	0.24	ug/l	
95-94-3	1,2,4,5-Tetrachlorobenzene	ND	2.2	0.41	ug/l	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
367-12-4	2-Fluorophenol	60%		14-88%
4165-62-2	Phenol-d5	39%		10-110%

ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound



Report of Analysis

Client Sample ID:	RA-18 GWS	Date Sampled:	05/02/16
Lab Sample ID:	JC19473-2	Date Received:	05/03/16
Matrix:	AQ - Ground Water	Percent Solids:	n/a
Method:	SW846 8270D SW846 3510C		
Project:	BMSMC, Building 5 Area, PR		

ABN TCL Special List

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
118-79-6	2,4,6-Tribromophenol	113%		39-149%
4165-60-0	Nitrobenzene-d5	111%		32-128%
321-60-8	2-Fluorobiphenyl	97%		35-119%
1718-51-0	Terphenyl-d14	93%		10-126%



ND = Not detected MDL = Method Detection Limit
RL = Reporting Limit
E = Indicates value exceeds calibration range

J = Indicates an estimated value
B = Indicates analyte found in associated method blank
N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID:	RA-18 GWS	Date Sampled:	05/02/16
Lab Sample ID:	JC19473-2	Date Received:	05/03/16
Matrix:	AQ - Ground Water	Percent Solids:	n/a
Method:	SW846 8270D BY SIM SW846 3510C		
Project:	BMSMC, Building 5 Area, PR		

	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	3M61048.D	1	05/04/16	LK	05/03/16	OP93555A	E3M2870
Run #2							

	Initial Volume	Final Volume
Run #1	900 ml	1.0 ml
Run #2		

CAS No.	Compound	Result	RL	MDL	Units	Q
91-20-3	Naphthalene	ND	0.11	0.033	ug/l	
123-91-1	1,4-Dioxane	2.62	0.11	0.054	ug/l	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
4165-60-0	Nitrobenzene-d5	88%		24-125%
321-60-8	2-Fluorobiphenyl	100%		19-127%
1718-51-0	Terphenyl-d14	106%		10-119%



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID: RA-18 GWS
 Lab Sample ID: JC19473-2
 Matrix: AQ - Ground Water
 Method: SW846 8081B SW846 3510C
 Project: BMSMC, Building 5 Area, PR

Date Sampled: 05/02/16
 Date Received: 05/03/16
 Percent Solids: n/a

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	6G34644.D	1	05/03/16	DS	05/03/16	OP93571	G6G1001
Run #2							

Run #	Initial Volume	Final Volume
Run #1	280 ml	2.0 ml
Run #2		

Pesticide TCL List

CAS No.	Compound	Result	RL	MDL	Units	Q
309-00-2	Aldrin	ND	0.0071	0.0043	ug/l	
319-84-6	alpha-BHC	ND	0.0071	0.0043	ug/l	
319-85-7	beta-BHC	ND	0.0071	0.0041	ug/l	
319-86-8	delta-BHC	ND	0.0071	0.0033	ug/l	
58-89-9	gamma-BHC (Lindane)	ND	0.0071	0.0020	ug/l	
5103-71-9	alpha-Chlordane	ND	0.0071	0.0033	ug/l	
5103-74-2	gamma-Chlordane	ND	0.0071	0.0033	ug/l	
60-57-1	Dieldrin	ND	0.0071	0.0026	ug/l	
72-54-8	4,4'-DDD	ND	0.0071	0.0027	ug/l	
72-55-9	4,4'-DDE	ND	0.0071	0.0044	ug/l	
50-29-3	4,4'-DDT	ND	0.0071	0.0035	ug/l	
72-20-8	Endrin	ND	0.0071	0.0036	ug/l	
1031-07-8	Endosulfan sulfate	ND	0.0071	0.0038	ug/l	
7421-93-4	Endrin aldehyde	ND	0.0071	0.0037	ug/l	
53494-70-5	Endrin ketone	ND	0.0071	0.0036	ug/l	
959-98-8	Endosulfan-I	ND	0.0071	0.0035	ug/l	
33213-65-9	Endosulfan-II	ND	0.0071	0.0031	ug/l	
76-44-8	Heptachlor	ND	0.0071	0.0027	ug/l	
1024-57-3	Heptachlor epoxide	ND	0.0071	0.0047	ug/l	
72-43-5	Methoxychlor	ND	0.014	0.0041	ug/l	
8001-35-2	Toxaphene	ND	0.18	0.13	ug/l	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
877-09-8	Tetrachloro-m-xylene	77%		26-132%
877-09-8	Tetrachloro-m-xylene	78%		26-132%
2051-24-3	Decachlorobiphenyl	103%		10-118%
2051-24-3	Decachlorobiphenyl	84%		10-118%



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID: RA-18 GWS		Date Sampled: 05/02/16	
Lab Sample ID: JC19473-2A		Date Received: 05/03/16	
Matrix: AQ - Ground Water		Percent Solids: n/a	
Method: SW846 8015C			
Project: BMSMC, Building 5 Area, PR			

Run #1 ^a	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #2	XY064115.D	1	05/05/16	AFL	n/a	n/a	F:GXY2774

CAS No.	Compound	Result	RL	MDL	Units	Q
64-17-5	Ethanol	ND	5.0	1.0	mg/l	
78-83-1	Isobutyl Alcohol	ND	5.0	1.0	mg/l	
67-63-0	Isopropyl Alcohol	ND	5.0	1.0	mg/l	
71-23-8	n-Propyl Alcohol	ND	5.0	1.0	mg/l	
71-36-3	n-Butyl Alcohol	ND	5.0	1.0	mg/l	
67-56-1	Methanol	ND	5.0	1.0	mg/l	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
111-27-3	Hexanol	98%		73-123%

(a) Analysis performed at Accutest Laboratories, Orlando FL.



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 3

Client Sample ID:	RA-18 (12.5-13.5)	Date Sampled:	05/02/16
Lab Sample ID:	JC19473-3	Date Received:	05/03/16
Matrix:	SO - Soil	Percent Solids:	76.1
Method:	SW846 8270D SW846 3546		
Project:	BMSMC, Building 5 Area, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	3E83106.D	1	05/11/16	AN	05/04/16	OP93597	E3E3643
Run #2							

Run #	Initial Weight	Final Volume
Run #1	30.5 g	1.0 ml
Run #2		

AEN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units	Q
95-57-8	2-Chlorophenol	ND	86	30	ug/kg	
59-50-7	4-Chloro-3-methyl phenol	ND	220	39	ug/kg	
120-83-2	2,4-Dichlorophenol	ND	220	42	ug/kg	
105-67-9	2,4-Dimethylphenol	ND	220	100	ug/kg	
51-28-5	2,4-Dinitrophenol	ND	220	170	ug/kg	
534-52-1	4,6-Dinitro-o-cresol	ND	220	63	ug/kg	
95-48-7	2-Methylphenol	ND	86	49	ug/kg	
	3&4-Methylphenol	ND	86	49	ug/kg	
88-75-5	2-Nitrophenol	ND	220	37	ug/kg	
100-02-7	4-Nitrophenol	ND	430	120	ug/kg	
87-86-5	Pentachlorophenol	ND	220	37	ug/kg	
108-95-2	Phenol	ND	86	34	ug/kg	
58-90-2	2,3,4,6-Tetrachlorophenol	ND	220	56	ug/kg	
95-95-4	2,4,5-Trichlorophenol	ND	220	62	ug/kg	
88-06-2	2,4,6-Trichlorophenol	ND	220	47	ug/kg	
83-32-9	Acenaphthene	ND	43	8.1	ug/kg	
208-96-8	Acenaphthylene	ND	43	5.9	ug/kg	
98-86-2	Acetophenone	ND	220	22	ug/kg	
120-12-7	Anthracene	ND	43	19	ug/kg	
1912-24-9	Atrazine	ND	86	14	ug/kg	
56-55-3	Benzo(a)anthracene	ND	43	6.6	ug/kg	
50-32-8	Benzo(a)pyrene	ND	43	11	ug/kg	
205-99-2	Benzo(b)fluoranthene	ND	43	9.0	ug/kg	
191-24-2	Benzo(g,h,i)perylene	ND	43	12	ug/kg	
207-08-9	Benzo(k)fluoranthene	ND	43	12	ug/kg	
101-55-3	4-Bromophenyl phenyl ether	ND	86	20	ug/kg	
85-68-7	Butyl benzyl phthalate	ND	86	25	ug/kg	
92-52-4	1,1'-Biphenyl	ND	86	12	ug/kg	
100-52-7	Benzaldehyde	ND	220	10	ug/kg	
91-58-7	2-Chloronaphthalene	ND	86	8.1	ug/kg	
106-47-8	4-Chloroaniline	ND	220	12	ug/kg	
86-74-8	Carbazole	ND	86	8.4	ug/kg	



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

Report of Analysis

Client Sample ID: RA-18 (12.5-13.5)
 Lab Sample ID: JC19473-3
 Matrix: SO - Soil
 Method: SW846 8270D SW846 3546
 Project: BSMC, Building 5 Area, PR

Date Sampled: 05/02/16
 Date Received: 05/03/16
 Percent Solids: 76.1

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units	Q
105-60-2	Caprolactam	ND	86	31	ug/kg	
218-01-9	Chrysene	ND	43	11	ug/kg	
111-91-1	bis(2-Chloroethoxy)methane	ND	86	9.1	ug/kg	
111-44-4	bis(2-Chloroethyl)ether	ND	86	25	ug/kg	
108-60-1	bis(2-Chloroisopropyl)ether	ND	86	17	ug/kg	
7005-72-3	4-Chlorophenyl phenyl ether	ND	86	11	ug/kg	
121-14-2	2,4-Dinitrotoluene	ND	43	16	ug/kg	
606-20-2	2,6-Dinitrotoluene	ND	43	13	ug/kg	
91-94-1	3,3'-Dichlorobenzidine	ND	86	58	ug/kg	
53-70-3	Dibenzo(a,h)anthracene	ND	43	16	ug/kg	
132-64-9	Dibenzofuran	ND	86	7.8	ug/kg	
84-74-2	Di-n-butyl phthalate	ND	86	27	ug/kg	
117-84-0	Di-n-octyl phthalate	ND	86	28	ug/kg	
84-66-2	Diethyl phthalate	ND	86	9.7	ug/kg	
131-11-3	Dimethyl phthalate	ND	86	8.8	ug/kg	
117-81-7	bis(2-Ethylhexyl)phthalate	ND	86	9.4	ug/kg	
206-44-0	Fluoranthene	ND	43	20	ug/kg	
86-73-7	Fluorene	ND	43	17	ug/kg	
118-74-1	Hexachlorobenzene	ND	86	11	ug/kg	
87-68-3	Hexachlorobutadiene	ND	43	17	ug/kg	
77-47-4	Hexachlorocyclopentadiene	ND	430	26	ug/kg	
67-72-1	Hexachloroethane	ND	220	33	ug/kg	
193-39-5	Indeno(1,2,3-cd)pyrene	ND	43	14	ug/kg	
78-59-1	Isophorone	ND	86	9.6	ug/kg	
90-12-0	1-Methylnaphthalene	ND	86	7.7	ug/kg	
91-57-6	2-Methylnaphthalene	ND	86	35	ug/kg	
88-74-4	2-Nitroaniline	ND	220	31	ug/kg	
99-09-2	3-Nitroaniline	ND	220	15	ug/kg	
100-01-6	4-Nitroaniline	ND	220	17	ug/kg	
98-95-3	Nitrobenzene	ND	86	19	ug/kg	
621-64-7	N-Nitroso-di-n-propylamine	ND	86	20	ug/kg	
86-30-6	N-Nitrosodiphenylamine	ND	220	25	ug/kg	
85-01-8	Phenanthrene	ND	43	10	ug/kg	
129-00-0	Pyrene	ND	43	7.5	ug/kg	
95-94-3	1,2,4,5-Tetrachlorobenzene	ND	220	9.7	ug/kg	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
367-12-4	2-Fluorophenol	48%		30-106%
4165-62-2	Phenol-d5	47%		30-106%

ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound



Report of Analysis

Client Sample ID: RA-18 (12.5-13.5)
Lab Sample ID: JC19473-3
Matrix: SO - Soil
Method: SW846 8270D SW846 3546
Project: BSMC, Building 5 Area, PR

Date Sampled: 05/02/16
Date Received: 05/03/16
Percent Solids: 76.1

ABN TCL Special List

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
118-79-6	2,4,6-Tribromophenol	54%		24-140%
4165-60-0	Nitrobenzene-d5	50%		26-122%
321-60-8	2-Fluorobiphenyl	59%		36-112%
1718-51-0	Terphenyl-d14	61%		36-132%



ND = Not detected MDL = Method Detection Limit
RL = Reporting Limit
E = Indicates value exceeds calibration range

J = Indicates an estimated value
B = Indicates analyte found in associated method blank
N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID:	RA-18 (12.5-13.5)	Date Sampled:	05/02/16
Lab Sample ID:	JC19473-3	Date Received:	05/03/16
Matrix:	SO - Soil	Percent Solids:	76.1
Method:	SW846 8270D BY SIM SW846 3546		
Project:	BMSMC, Building 5 Area, PR		

	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	3M61096.D	1	05/05/16	LK	05/04/16	OP93597A	E3M2872
Run #2							

	Initial Weight	Final Volume
Run #1	30.5 g	1.0 ml
Run #2		

CAS No.	Compound	Result	RL	MDL	Units	Q
123-91-1	1,4-Dioxane	ND	4.3	0.87	ug/kg	
91-20-3	Naphthalene	ND	4.3	0.53	ug/kg	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
4165-60-0	Nitrobenzene-d5	57%		15-138%
321-60-8	2-Fluorobiphenyl	69%		12-148%
1718-51-0	Terphenyl-d14	74%		10-157%



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID:	RA-18 (12.5-13.5)	Date Sampled:	05/02/16
Lab Sample ID:	JC19473-3	Date Received:	05/03/16
Matrix:	SO - Soil	Percent Solids:	76.1
Method:	SW846 8015C MOD		
Project:	BMSMC, Building 5 Area, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	XY064123.D	1	05/06/16	AFL	n/a	n/a	F:GXY2775
Run #2							

Run #	Initial Weight	Final Volume
Run #1	4.77 g	10.0 ml
Run #2		

CAS No.	Compound	Result	RL	MDL	Units	Q
64-17-5	Ethanol	ND	14	2.8	mg/kg	
78-83-1	Isobutyl Alcohol	ND	14	2.8	mg/kg	
67-63-0	Isopropyl Alcohol	ND	14	2.8	mg/kg	
71-23-8	n-Propyl Alcohol	ND	14	2.8	mg/kg	
71-36-3	n-Butyl Alcohol	ND	14	2.8	mg/kg	
67-56-1	Methanol	ND	14	2.8	mg/kg	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
111-27-3	Hexanol	130% ^b		69-121%

(a) Sample was received in a bulk container but was not preserved within 48 hours of sampling. Analysis performed at Accutest Laboratories, Orlando FL.

(b) Outside control limits; however, sample is ND.



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID:	RA-18 (12.5-13.5)	Date Sampled:	05/02/16
Lab Sample ID:	JC19473-3	Date Received:	05/03/16
Matrix:	SO - Soil	Percent Solids:	76.1
Method:	SW846 8081B SW846 3546		
Project:	BMSMC, Building 5 Area, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	4G68007.D	1	05/06/16	BP	05/06/16	OP93683	G4G1782
Run #2							

Run #	Initial Weight	Final Volume
Run #1	15.3 g	10.0 ml
Run #2		

Pesticide TCL List

CAS No.	Compound	Result	RL	MDL	Units	Q
309-00-2	Aldrin	ND	0.86	0.77	ug/kg	
319-84-6	alpha-BHC	ND	0.86	0.57	ug/kg	
319-85-7	beta-BHC	ND	0.86	0.53	ug/kg	
319-86-8	delta-BHC	ND	0.86	0.34	ug/kg	
58-89-9	gamma-BHC (Lindane)	ND	0.86	0.39	ug/kg	
5103-71-9	alpha-Chlordane	ND	0.86	0.46	ug/kg	
5103-74-2	gamma-Chlordane	ND	0.86	0.65	ug/kg	
60-57-1	Dieldrin	ND	0.86	0.67	ug/kg	
72-54-8	4,4'-DDD	ND	0.86	0.32	ug/kg	
72-55-9	4,4'-DDE	ND	0.86	0.29	ug/kg	
50-29-3	4,4'-DDT	ND	0.86	0.33	ug/kg	
72-20-8	Endrin	ND	0.86	0.30	ug/kg	
1031-07-8	Endosulfan sulfate	ND	0.86	0.49	ug/kg	
7421-93-4	Endrin aldehyde	ND	0.86	0.64	ug/kg	
959-98-8	Endosulfan-I	ND	0.86	0.28	ug/kg	
33213-65-9	Endosulfan-II	ND	0.86	0.81	ug/kg	
76-44-8	Heptachlor	ND	0.86	0.71	ug/kg	
1024-57-3	Heptachlor epoxide	ND	0.86	0.35	ug/kg	
72-43-5	Methoxychlor	ND	1.7	0.48	ug/kg	
53494-70-5	Endrin ketone	ND	0.86	0.45	ug/kg	
8001-35-2	Toxaphene	ND	21	15	ug/kg	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
877-09-8	Tetrachloro-m-xylene	110%		24-136%
877-09-8	Tetrachloro-m-xylene	107%		24-136%
2051-24-3	Decachlorobiphenyl	139%		10-153%
2051-24-3	Decachlorobiphenyl	108%		10-153%



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS

GW
W/EB
SD - Accutest
ACCUTEST - NJ

CHAIN OF CUSTODY

SGS Accutest - Dayton
2333 Route 130, Dayton, NJ 08810
TEL 732-329-0200 FAX 732-329-3490
www.accutest.com

PAGE 1 OF 1 5-0

801219536121

Form CCR-0000-0

SGS Accutest Chain #

SGS Accutest Job #

JC19473

Company Name Anderson Mulholland Assoc. Inc.		Project Name BMS Release Assessment		Matrix Codes	
Street Address 2700 Westchester		City Humacao PR		Matrix Codes	
City Purchase NY		Billing Information (If different from Report to)		Matrix Codes	
Project Contact Terry Taylor		Company Name		Matrix Codes	
Phone # 914-251-0400		Client Purchase Order #		Matrix Codes	
Fax #		City		Matrix Codes	
State		State		Matrix Codes	
Zip		Zip		Matrix Codes	
Project Manager N. Rivera, T. Taylor, D. Lindstrom		Analyst		Matrix Codes	
SGS Account Number		Collection		Matrix Codes	
Field ID / Point of Collection		McGraw-Hill #		Matrix Codes	
Date		Time		Matrix Codes	
Sampled by		Method		Matrix Codes	
# of bottles		# of bottles		Matrix Codes	
HCl		HCl		Matrix Codes	
HNO3		HNO3		Matrix Codes	
H2SO4		H2SO4		Matrix Codes	
H3PO4		H3PO4		Matrix Codes	
H4SiO4		H4SiO4		Matrix Codes	
HNO2		HNO2		Matrix Codes	
HNO		HNO		Matrix Codes	
H2O2		H2O2		Matrix Codes	
H2S		H2S		Matrix Codes	
H2SO3		H2SO3		Matrix Codes	
H2S2O8		H2S2O8		Matrix Codes	
H2S2O7		H2S2O7		Matrix Codes	
H2S2O5		H2S2O5		Matrix Codes	
H2S2O3		H2S2O3		Matrix Codes	
H2S2O2		H2S2O2		Matrix Codes	
H2S2O		H2S2O		Matrix Codes	
H2S		H2S		Matrix Codes	
H2		H2		Matrix Codes	
H		H		Matrix Codes	
O		O		Matrix Codes	
N		N		Matrix Codes	
C		C		Matrix Codes	
S		S		Matrix Codes	
P		P		Matrix Codes	
K		K		Matrix Codes	
Ca		Ca		Matrix Codes	
Mg		Mg		Matrix Codes	
Na		Na		Matrix Codes	
Fe		Fe		Matrix Codes	
Cu		Cu		Matrix Codes	
Zn		Zn		Matrix Codes	
Mn		Mn		Matrix Codes	
Co		Co		Matrix Codes	
Ni		Ni		Matrix Codes	
Cd		Cd		Matrix Codes	
Cr		Cr		Matrix Codes	
M		M		Matrix Codes	
B		B		Matrix Codes	
Al		Al		Matrix Codes	
Si		Si		Matrix Codes	
Cl		Cl		Matrix Codes	
Br		Br		Matrix Codes	
I		I		Matrix Codes	
F		F		Matrix Codes	
Li		Li		Matrix Codes	
Be		Be		Matrix Codes	
B		B		Matrix Codes	
C		C		Matrix Codes	
N		N		Matrix Codes	
O		O		Matrix Codes	
F		F		Matrix Codes	
Cl		Cl		Matrix Codes	
Br		Br		Matrix Codes	
I		I		Matrix Codes	
S		S		Matrix Codes	
Se		Se		Matrix Codes	
Te		Te		Matrix Codes	
Pb		Pb		Matrix Codes	
Cu		Cu		Matrix Codes	
Zn		Zn		Matrix Codes	
Mn		Mn		Matrix Codes	
Co		Co		Matrix Codes	
Ni		Ni		Matrix Codes	
Cd		Cd		Matrix Codes	
Cr		Cr		Matrix Codes	
M		M		Matrix Codes	
B		B		Matrix Codes	
Al		Al		Matrix Codes	
Si		Si		Matrix Codes	
Cl		Cl		Matrix Codes	
Br		Br		Matrix Codes	
I		I		Matrix Codes	
F		F		Matrix Codes	
Li		Li		Matrix Codes	
Be		Be		Matrix Codes	
B		B		Matrix Codes	
C		C		Matrix Codes	
N		N		Matrix Codes	
O		O		Matrix Codes	
F		F		Matrix Codes	
Cl		Cl		Matrix Codes	
Br		Br		Matrix Codes	
I		I		Matrix Codes	
S		S		Matrix Codes	
Se		Se		Matrix Codes	
Te		Te		Matrix Codes	
Pb		Pb		Matrix Codes	
Cu		Cu		Matrix Codes	
Zn		Zn		Matrix Codes	
Mn		Mn		Matrix Codes	
Co		Co		Matrix Codes	
Ni		Ni		Matrix Codes	
Cd		Cd		Matrix Codes	
Cr		Cr		Matrix Codes	
M		M		Matrix Codes	
B		B		Matrix Codes	
Al		Al		Matrix Codes	
Si		Si		Matrix Codes	
Cl		Cl		Matrix Codes	
Br		Br		Matrix Codes	
I		I		Matrix Codes	
F		F		Matrix Codes	
Li		Li		Matrix Codes	
Be		Be		Matrix Codes	
B		B		Matrix Codes	
C		C		Matrix Codes	
N		N		Matrix Codes	
O		O		Matrix Codes	
F		F		Matrix Codes	
Cl		Cl		Matrix Codes	
Br		Br		Matrix Codes	
I		I		Matrix Codes	
S		S		Matrix Codes	
Se		Se		Matrix Codes	
Te		Te		Matrix Codes	
Pb		Pb		Matrix Codes	
Cu		Cu		Matrix Codes	
Zn		Zn		Matrix Codes	
Mn		Mn		Matrix Codes	
Co		Co		Matrix Codes	
Ni		Ni		Matrix Codes	
Cd		Cd		Matrix Codes	
Cr		Cr		Matrix Codes	
M		M		Matrix Codes	
B		B		Matrix Codes	
Al		Al		Matrix Codes	
Si		Si		Matrix Codes	
Cl		Cl		Matrix Codes	
Br		Br		Matrix Codes	
I		I		Matrix Codes	
F		F		Matrix Codes	
Li		Li		Matrix Codes	
Be		Be		Matrix Codes	
B		B		Matrix Codes	
C		C		Matrix Codes	
N		N		Matrix Codes	
O		O		Matrix Codes	
F		F		Matrix Codes	
Cl		Cl		Matrix Codes	
Br		Br		Matrix Codes	
I		I		Matrix Codes	
S		S		Matrix Codes	
Se		Se		Matrix Codes	
Te		Te		Matrix Codes	
Pb		Pb		Matrix Codes	
Cu		Cu		Matrix Codes	
Zn		Zn		Matrix Codes	
Mn		Mn		Matrix Codes	
Co		Co		Matrix Codes	
Ni		Ni		Matrix Codes	
Cd		Cd		Matrix Codes	
Cr		Cr		Matrix Codes	
M		M		Matrix Codes	
B		B		Matrix Codes	
Al		Al		Matrix Codes	
Si		Si		Matrix Codes	
Cl		Cl		Matrix Codes	
Br		Br		Matrix Codes	
I		I		Matrix Codes	
F		F		Matrix Codes	
Li		Li		Matrix Codes	
Be		Be		Matrix Codes	
B		B		Matrix Codes	
C		C		Matrix Codes	
N		N		Matrix Codes	
O		O		Matrix Codes	
F		F		Matrix Codes	
Cl		Cl		Matrix Codes	
Br		Br		Matrix Codes	
I		I		Matrix Codes	
S		S		Matrix Codes	
Se		Se		Matrix Codes	
Te		Te		Matrix Codes	
Pb		Pb		Matrix Codes	
Cu		Cu		Matrix Codes	
Zn		Zn		Matrix Codes	
Mn		Mn		Matrix Codes	
Co		Co		Matrix Codes	
Ni		Ni		Matrix Codes	
Cd		Cd		Matrix Codes	
Cr		Cr		Matrix Codes	
M		M		Matrix Codes	
B		B		Matrix Codes	
Al		Al		Matrix Codes	
Si		Si		Matrix Codes	
Cl		Cl		Matrix Codes	
Br		Br		Matrix Codes	
I		I		Matrix Codes	
F		F		Matrix Codes	
Li		Li		Matrix Codes	
Be		Be		Matrix Codes	
B		B		Matrix Codes	
C		C		Matrix Codes	
N		N		Matrix Codes	
O		O		Matrix Codes	
F		F		Matrix Codes	
Cl		Cl		Matrix Codes	
Br		Br		Matrix Codes	
I		I		Matrix Codes	
S		S		Matrix Codes	
Se		Se		Matrix Codes	
Te		Te		Matrix Codes	
Pb		Pb		Matrix Codes	
Cu		Cu		Matrix Codes	
Zn		Zn		Matrix Codes	
Mn		Mn		Matrix Codes	
Co		Co		Matrix Codes	
Ni		Ni		Matrix Codes	
Cd		Cd		Matrix Codes	
Cr		Cr		Matrix Codes	
M		M		Matrix Codes	
B		B		Matrix Codes	
Al		Al		Matrix Codes	
Si		Si		Matrix Codes	
Cl		Cl		Matrix Codes	
Br		Br		Matrix Codes	
I		I		Matrix Codes	
F		F		Matrix Codes	
Li		Li		Matrix Codes	
Be		Be		Matrix Codes	
B		B		Matrix Codes	
C		C		Matrix Codes	
N		N		Matrix Codes	
O		O		Matrix Codes	
F		F		Matrix Codes	
Cl		Cl		Matrix Codes	
Br		Br		Matrix Codes	
I		I		Matrix Codes	
S		S		Matrix Codes	
Se		Se		Matrix Codes	
Te		Te		Matrix Codes	
Pb		Pb		Matrix Codes	
Cu		Cu		Matrix Codes	
Zn		Zn		Matrix Codes	
Mn		Mn		Matrix Codes	
Co		Co		Matrix Codes	
Ni		Ni		Matrix Codes	
Cd		Cd		Matrix Codes	
Cr		Cr		Matrix Codes	
M		M		Matrix Codes	
B		B		Matrix Codes	
Al		Al		Matrix Codes	
Si		Si		Matrix Codes	
Cl		Cl		Matrix Codes	
Br		Br		Matrix Codes	
I		I		Matrix Codes	
F		F		Matrix Codes	
Li		Li		Matrix Codes	
Be		Be		Matrix Codes	
B		B		Matrix Codes	
C		C		Matrix Codes	
N		N		Matrix Codes	
O		O		Matrix Codes	
F		F		Matrix Codes	
Cl		Cl		Matrix Codes	
Br		Br		Matrix Codes	
I		I		Matrix Codes	
S		S		Matrix Codes	
Se		Se		Matrix Codes	
Te		Te		Matrix Codes	
Pb		Pb		Matrix Codes	
Cu		Cu		Matrix Codes	
Zn		Zn		Matrix Codes	
Mn		Mn		Matrix Codes	
Co		Co		Matrix Codes	
Ni		Ni		Matrix Codes	
Cd		Cd		Matrix Codes	
Cr		Cr		Matrix Codes	
M		M		Matrix Codes	
B		B		Matrix Codes	
Al		Al		Matrix Codes	
Si		Si		Matrix Codes	
Cl		Cl		Matrix Codes	
Br		Br		Matrix Codes	
I		I		Matrix Codes	
F		F		Matrix Codes	
Li		Li		Matrix Codes	
Be		Be		Matrix Codes	
B		B		Matrix Codes	
C		C		Matrix Codes	
N		N		Matrix Codes	
O		O		Matrix Codes	
F		F		Matrix Codes	
Cl		Cl		Matrix Codes	
Br		Br		Matrix Codes	
I		I		Matrix Codes	
S		S		Matrix Codes	
Se		Se		Matrix Codes	
Te		Te		Matrix Codes	
Pb		P			

EXECUTIVE NARRATIVE

SDG No: **JC19473** Laboratory: **Accutest, New Jersey**
Analysis: **SW846-8270D** Number of Samples: **3**
Location: **BMSMC, Building 5 Area**
Humacao, PR

SUMMARY: Three (3) samples were analyzed for the ABN TCL list following method SW846-8270D; Naphthalene and 1,4-Dioxane were also analyzed by SW846-8270D using the selective ion monitoring (SIM) technique. The sample results were assessed according to USEPA data validation guidance documents in the following order of precedence: EPA Hazardous Waste Support Section, SOP HW-35A, July 2015 –Revision 0. *Semivolatile Data Validation*. The QC criteria and data validation actions listed on the data review worksheets are from the primary guidance document, unless otherwise noted.

Results are valid and can be used for decision making purposes.

Critical issues: **None**
Major: **None**
Minor: **None**

Critical findings: **None**
Major findings: **None**


Minor findings:

1. Initial and continuing calibration verifications meet the required criteria. Analytes not meeting the method % difference criteria meet the guidance document performance criteria for continuing calibration verification of ± 25 or 40 %, no action taken. No closing calibration verification included in data package. No action taken, professional judgment.
2. Analytes not meeting the continuing calibration verification criteria of the guidance document were qualified UJ in samples JC19473-1; JC19473-2; and JC19473-3.
3. MS/MSD % recoveries outside control limits for several analytes. No action taken, MS/MSD results apply to the unspike sample. QC samples were from another project.

COMMENTS: Results are valid and can be used for decision making purposes.

Reviewers Name: **Rafael Infante**
Chemist License 1888

Signature:

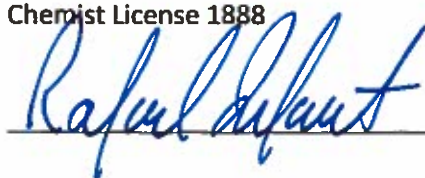

Date: **May 22, 2016**

SAMPLE ORGANIC DATA SAMPLE SUMMARY

Sample ID: JC19473-1

Sample location: BMSMC Building 5 Area

Sampling date: 4/30/2016

Matrix: AQ - Equipment Blank

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
2-Chlorophenol	5.6	ug/l	1	-	U	Yes
4-Chloro-3-methyl phenol	5.6	ug/l	1	-	U	Yes
2,4-Dichlorophenol	2.2	ug/l	1	-	U	Yes
2,4-Dimethylphenol	5.6	ug/l	1	-	U	Yes
2,4-Dinitrophenol	11	ug/l	1	-	U	Yes
4,6-Dinitro-o-cresol	5.6	ug/l	1	-	U	Yes
2-Methylphenol	2.2	ug/l	1	-	U	Yes
3&4-Methylphenol	2.2	ug/l	1	-	U	Yes
2-Nitrophenol	5.6	ug/l	1	-	U	Yes
4-Nitrophenol	11	ug/l	1	-	U	Yes
Pentachlorophenol	5.6	ug/l	1	-	U	Yes
Phenol	2.2	ug/l	1	-	U	Yes
2,3,4,6-Tetrachlorophenol	5.6	ug/l	1	-	U	Yes
2,4,5-Trichlorophenol	5.6	ug/l	1	-	U	Yes
2,4,6-Trichlorophenol	5.6	ug/l	1	-	U	Yes
Acenaphthene	1.1	ug/l	1	-	U	Yes
Acenaphthylene	1.1	ug/l	1	-	U	Yes
Acetophenone	2.2	ug/l	1	-	U	Yes
Anthracene	1.1	ug/l	1	-	U	Yes
Atrazine	2.2	ug/l	1	-	U	Yes
Benzaldehyde	5.6	ug/l	1	-	U	Yes
Benzo(a)anthracene	1.1	ug/l	1	-	U	Yes
Benzo(a)pyrene	1.1	ug/l	1	-	U	Yes
Benzo(b)fluoranthene	1.1	ug/l	1	-	U	Yes
Benzo(g,h,i)perylene	1.1	ug/l	1	-	U	Yes
Benzo(k)fluoranthene	1.1	ug/l	1	-	U	Yes
4-Bromophenyl phenyl ether	2.2	ug/l	1	-	U	Yes
Butyl benzyl phthalate	2.2	ug/l	1	-	U	Yes
1,1'-Biphenyl	1.1	ug/l	1	-	U	Yes
2-Chloronaphthalene	2.2	ug/l	1	-	U	Yes
4-Chloroaniline	5.6	ug/l	1	-	U	Yes
Carbazole	1.1	ug/l	1	-	U	Yes
Caprolactam	2.2	ug/l	1	-	U	Yes
Chrysene	1.1	ug/l	1	-	U	Yes
bis(2-Chloroethoxy)methane	2.2	ug/l	1	-	U	Yes
bis(2-Chloroethyl)ether	2.2	ug/l	1	-	U	Yes

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
bis(2-Chloroisopropyl)ether	2.2	ug/l	1	-	U	Yes
4-Chlorophenyl phenyl ether	2.2	ug/l	1	-	U	Yes
2,4-Dinitrotoluene	1.1	ug/l	1	-	U	Yes
2,6-Dinitrotoluene	1.1	ug/l	1	-	U	Yes
3,3'-Dichlorobenzidine	2.2	ug/l	1	-	U	Yes
Dibenzo(a,h)anthracene	1.1	ug/l	1	-	U	Yes
Dibenzofuran	5.6	ug/l	1	-	U	Yes
Di-n-butyl phthalate	2.2	ug/l	1	-	U	Yes
Di-n-octyl phthalate	2.2	ug/l	1	-	U	Yes
Diethyl phthalate	2.2	ug/l	1	-	U	Yes
Dimethyl phthalate	2.2	ug/l	1	-	U	Yes
bis(2-Ethylhexyl)phthalate	2.2	ug/l	1	-	U	Yes
Fluoranthene	1.1	ug/l	1	-	U	Yes
Fluorene	1.1	ug/l	1	-	U	Yes
Hexachlorobenzene	1.1	ug/l	1	-	U	Yes
Hexachlorobutadiene	1.1	ug/l	1	-	U	Yes
Hexachlorocyclopentadiene	11	ug/l	1	-	U	Yes
Hexachloroethane	2.2	ug/l	1	-	U	Yes
Indeno(1,2,3-cd)pyrene	1.1	ug/l	1	-	U	Yes
Isophorone	2.2	ug/l	1	-	U	Yes
1-Methylnaphthalene	1.1	ug/l	1	-	U	Yes
2-Methylnaphthalene	1.1	ug/l	1	-	U	Yes
2-Nitroaniline	5.6	ug/l	1	-	UJ	Yes
3-Nitroaniline	5.6	ug/l	1	-	U	Yes
4-Nitroaniline	5.6	ug/l	1	-	U	Yes
Nitrobenzene	2.2	ug/l	1	-	U	Yes
N-Nitroso-di-n-propylamine	2.2	ug/l	1	-	UJ	Yes
Nitrosodiphenylamine	5.6	ug/l	1	-	U	Yes
Phenanthrene	1.1	ug/l	1	-	U	Yes
Pyrene	1.1	ug/l	1	-	U	Yes
1,2,4,5-Tetrachlorobenzene	2.2	ug/l	1	-	UJ	Yes

METHOD: 8270D (SIM)

Naphthalene	0.11	ug/l	1	-	U	Yes
1,4-Dioxane	0.11	ug/l	1	-	U	Yes

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
--------------	--------	-------	-----------------	----------	------------	------------

Sample ID: JC19473-2

Sample location: BSMC Building 5 Area

Sampling date: 5/3/2016

Matrix: Groundwater

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
2-Chlorophenol	5.6	ug/l	1	-	U	Yes
4-Chloro-3-methyl phenol	5.6	ug/l	1	-	U	Yes
2,4-Dichlorophenol	2.2	ug/l	1	-	U	Yes
2,4-Dimethylphenol	5.6	ug/l	1	-	U	Yes
2,4-Dinitrophenol	11	ug/l	1	-	U	Yes
4,6-Dinitro-o-cresol	5.6	ug/l	1	-	U	Yes
2-Methylphenol	2.2	ug/l	1	-	U	Yes
3&4-Methylphenol	2.2	ug/l	1	-	U	Yes
2-Nitrophenol	5.6	ug/l	1	-	U	Yes
4-Nitrophenol	11	ug/l	1	-	U	Yes
Pentachlorophenol	5.6	ug/l	1	-	U	Yes
Phenol	2.2	ug/l	1	-	U	Yes
2,3,4,6-Tetrachlorophenol	5.6	ug/l	1	-	U	Yes
2,4,5-Trichlorophenol	5.6	ug/l	1	-	U	Yes
2,4,6-Trichlorophenol	5.6	ug/l	1	-	U	Yes
Acenaphthene	1.1	ug/l	1	-	U	Yes
Acenaphthylene	1.1	ug/l	1	-	U	Yes
Acetophenone	2.2	ug/l	1	-	U	Yes
Anthracene	1.1	ug/l	1	-	U	Yes
Atrazine	2.2	ug/l	1	-	U	Yes
Benzaldehyde	5.6	ug/l	1	-	U	Yes
Benzo(a)anthracene	1.1	ug/l	1	-	U	Yes
Benzo(a)pyrene	1.1	ug/l	1	-	U	Yes
Benzo(b)fluoranthene	1.1	ug/l	1	-	U	Yes
Benzo(g,h,i)perylene	1.1	ug/l	1	-	U	Yes
Benzo(k)fluoranthene	1.1	ug/l	1	-	U	Yes
4-Bromophenyl phenyl ether	2.2	ug/l	1	-	U	Yes
Butyl benzyl phthalate	2.2	ug/l	1	-	U	Yes
1,1'-Biphenyl	1.1	ug/l	1	-	U	Yes
2-Chloronaphthalene	2.2	ug/l	1	-	U	Yes
4-Chloroaniline	5.6	ug/l	1	-	U	Yes
Carbazole	1.1	ug/l	1	-	U	Yes
Caprolactam	2.2	ug/l	1	-	U	Yes
Chrysene	1.1	ug/l	1	-	U	Yes
bis(2-Chloroethoxy)methane	2.2	ug/l	1	-	U	Yes
bis(2-Chloroethyl)ether	2.2	ug/l	1	-	U	Yes

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
bis(2-Chloroisopropyl)ether	2.2	ug/l	1	-	U	Yes
4-Chlorophenyl phenyl ether	2.2	ug/l	1	-	U	Yes
2,4-Dinitrotoluene	1.1	ug/l	1	-	U	Yes
2,6-Dinitrotoluene	1.1	ug/l	1	-	U	Yes
3,3'-Dichlorobenzidine	2.2	ug/l	1	-	U	Yes
Dibenzo(a,h)anthracene	1.1	ug/l	1	-	U	Yes
Dibenzofuran	5.6	ug/l	1	-	U	Yes
Di-n-butyl phthalate	2.2	ug/l	1	-	U	Yes
Di-n-octyl phthalate	2.2	ug/l	1	-	U	Yes
Diethyl phthalate	2.2	ug/l	1	-	U	Yes
Dimethyl phthalate	2.2	ug/l	1	-	U	Yes
bis(2-Ethylhexyl)phthalate	2.2	ug/l	1	-	U	Yes
Fluoranthene	1.1	ug/l	1	-	U	Yes
Fluorene	1.1	ug/l	1	-	U	Yes
Hexachlorobenzene	1.1	ug/l	1	-	U	Yes
Hexachlorobutadiene	1.1	ug/l	1	-	U	Yes
Hexachlorocyclopentadiene	11	ug/l	1	-	U	Yes
Hexachloroethane	2.2	ug/l	1	-	U	Yes
Indeno(1,2,3-cd)pyrene	1.1	ug/l	1	-	U	Yes
Isophorone	2.2	ug/l	1	-	U	Yes
1-Methylnaphthalene	1.1	ug/l	1	-	U	Yes
2-Methylnaphthalene	1.1	ug/l	1	-	U	Yes
2-Nitroaniline	5.6	ug/l	1	-	UJ	Yes
3-Nitroaniline	5.6	ug/l	1	-	U	Yes
4-Nitroaniline	5.6	ug/l	1	-	U	Yes
Nitrobenzene	2.2	ug/l	1	-	U	Yes
N-Nitroso-di-n-propylamine	2.2	ug/l	1	-	UJ	Yes
Nitrosodiphenylamine	5.6	ug/l	1	-	U	Yes
Phenanthrene	1.1	ug/l	1	-	U	Yes
Pyrene	1.1	ug/l	1	-	U	Yes
1,2,4,5-Tetrachlorobenzene	2.2	ug/l	1	-	UJ	Yes

METHOD: 8270D (SIM)

Naphthalene	0.11	ug/l	1	-	U	Yes
1,4-Dioxane	2.62	ug/l	100	-	-	Yes

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
--------------	--------	-------	-----------------	----------	------------	------------

Sample ID: JC19023-3

Sample location: BMSMC Building 5 Area

Sampling date: 5/2/2016

Matrix: Soil

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
2-Chlorophenol	86	ug/kg	1	-	U	Yes
4-Chloro-3-methyl phenol	220	ug/kg	1	-	U	Yes
2,4-Dichlorophenol	220	ug/kg	1	-	U	Yes
2,4-Dimethylphenol	220	ug/kg	1	-	U	Yes
2,4-Dinitrophenol	220	ug/kg	1	-	U	Yes
4,6-Dinitro-o-cresol	220	ug/kg	1	-	U	Yes
2-Methylphenol	86	ug/kg	1	-	U	Yes
3&4-Methylphenol	86	ug/kg	1	-	U	Yes
2-Nitrophenol	220	ug/kg	1	-	U	Yes
4-Nitrophenol	10	ug/kg	1	-	U	Yes
Pentachlorophenol	220	ug/kg	1	-	U	Yes
Phenol	86	ug/kg	1	-	U	Yes
2,3,4,6-Tetrachlorophenol	220	ug/kg	1	-	U	Yes
2,4,5-Trichlorophenol	220	ug/kg	1	-	U	Yes
2,4,6-Trichlorophenol	220	ug/kg	1	-	U	Yes
Acenaphthene	43	ug/kg	1	-	U	Yes
Acenaphthylene	43	ug/kg	1	-	U	Yes
Acetophenone	220	ug/kg	1	-	U	Yes
Anthracene	43	ug/kg	1	-	U	Yes
Atrazine	86	ug/kg	1	-	U	Yes
Benzo(a)anthracene	43	ug/kg	1	-	U	Yes
Benzo(a)pyrene	43	ug/kg	1	-	U	Yes
Benzo(b)fluoranthene	43	ug/kg	1	-	U	Yes
Benzo(g,h,i)perylene	43	ug/kg	1	-	U	Yes
Benzo(k)fluoranthene	43	ug/kg	1	-	U	Yes
4-Bromophenyl phenyl ether	86	ug/kg	1	-	U	Yes
Butyl benzyl phthalate	86	ug/kg	1	-	U	Yes
1,1'-Biphenyl	86	ug/kg	1	-	U	Yes
Benzaldehyde	220	ug/kg	1	-	U	Yes
2-Chloronaphthalene	86	ug/kg	1	-	U	Yes
4-Chloroaniline	220	ug/kg	1	-	U	Yes
Carbazole	86	ug/kg	1	-	U	Yes
Caprolactam	86	ug/kg	1	-	U	Yes
Chrysene	43	ug/kg	1	-	U	Yes
bis(2-Chloroethoxy)methane	86	ug/kg	1	-	U	Yes
bis(2-Chloroethyl)ether	86	ug/kg	1	-	U	Yes

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
bis(2-Chloroisopropyl)ether	86	ug/kg	1	-	U	Yes
4-Chlorophenyl phenyl ether	86	ug/kg	1	-	U	Yes
2,4-Dinitrotoluene	43	ug/kg	1	-	U	Yes
2,6-Dinitrotoluene	43	ug/kg	1	-	U	Yes
3,3'-Dichlorobenzidine	86	ug/kg	1	-	U	Yes
Dibenzo(a,h)anthracene	43	ug/kg	1	-	U	Yes
Dibenzofuran	86	ug/kg	1	-	U	Yes
Di-n-butyl phthalate	86	ug/kg	1	-	U	Yes
Di-n-octyl phthalate	86	ug/kg	1	-	U	Yes
Diethyl phthalate	86	ug/kg	1	-	U	Yes
Dimethyl phthalate	86	ug/kg	1	-	U	Yes
bis(2-Ethylhexyl)phthalate	86	ug/kg	1	-	U	Yes
Fluoranthene	43	ug/kg	1	-	U	Yes
Fluorene	43	ug/kg	1	-	U	Yes
Hexachlorobenzene	86	ug/kg	1	-	U	Yes
Hexachlorobutadiene	43	ug/kg	1	-	U	Yes
Hexachlorocyclopentadiene	430	ug/kg	1	-	U	Yes
Hexachloroethane	220	ug/kg	1	-	U	Yes
Indeno(1,2,3-cd)pyrene	43	ug/kg	1	-	U	Yes
Isophorone	86	ug/kg	1	-	U	Yes
1-Methylnaphthalene	86	ug/kg	1	-	U	Yes
2-Methylnaphthalene	86	ug/kg	1	-	U	Yes
2-Nitroaniline	220	ug/kg	1	-	U	Yes
3-Nitroaniline	220	ug/kg	1	-	U	Yes
4-Nitroaniline	220	ug/kg	1	-	U	Yes
Nitrobenzene	86	ug/kg	1	-	U	Yes
N-Nitroso-di-n-propylamine	86	ug/kg	1	-	U	Yes
N-Nitrosodiphenylamine	220	ug/kg	1	-	U	Yes
Phenanthrene	43	ug/kg	1	-	U	Yes
Pyrene	43	ug/kg	1	-	U	Yes
1,2,4,5-Tetrachlorobenzene	220	ug/kg	1	-	UJ	Yes

METHOD: 8270D (SIM)

Naphthalene	4.3	ug/kg	1	-	J	Yes
1,4-Dioxane	4.3	ug/kg	1	-	R	Yes

DATA REVIEW WORKSHEETS

Project Number: JC19473
 Date: April 30-May 02, 2016
 Shipping Date: May 02, 2016
 EPA Region: 2

REVIEW OF SEMIVOLATILE ORGANIC PACKAGE

The following guidelines for evaluating volatile organics were created to delineate required validation actions. This document will assist the reviewer in using professional judgment to make more informed decision and in better serving the needs of the data users. The sample results were assessed according to USEPA data validation guidance documents in the following order of precedence: EPA Hazardous Waste Support Section, SOP HW-35A, July 2015 –Revision 0. *Semivolatile Data Validation*. The QC criteria and data validation actions listed on the data review worksheets are from the primary guidance document, unless otherwise noted.

The hardcopied (laboratory name) Accutest data package received has been reviewed and the quality control and performance data summarized. The data review for SVOCs included:

Lab. Project/SDG No.: JC19473 Sample matrix: Soil/Groundwater
 No. of Samples: 3_Full_scan/3_SIM

Trip blank No.: -

Field blank No.: -

Equipment blank No.: JC19473-1

Field duplicate No.: -

<input checked="" type="checkbox"/> Data Completeness	<input checked="" type="checkbox"/> Laboratory Control Spikes
<input checked="" type="checkbox"/> Holding Times	<input checked="" type="checkbox"/> Field Duplicates
<input checked="" type="checkbox"/> GC/MS Tuning	<input checked="" type="checkbox"/> Calibrations
<input checked="" type="checkbox"/> Internal Standard Performance	<input checked="" type="checkbox"/> Compound Identifications
<input checked="" type="checkbox"/> Blanks	<input checked="" type="checkbox"/> Compound Quantitation
<input checked="" type="checkbox"/> Surrogate Recoveries	<input checked="" type="checkbox"/> Quantitation Limits
<input checked="" type="checkbox"/> Matrix Spike/Matrix Spike Duplicate	

Overall Comments: ABN_TCL_list_by_method_SW846-8270D;_Naphthalene_and_1,4-Dioxane_
_analyzed_by_method_SW846-8270D_(SIM)

Definition of Qualifiers:

J- Estimated results
 U- Compound not detected
 R- Rejected data
 UJ- Estimated nondetect

Reviewer: 
 Date: May 22, 2016

DATA COMPLETENESS

DATE RECEIVED

This image shows a single sheet of white paper with horizontal blue or grey ruling lines. The lines are evenly spaced and run across the width of the page. There is no handwriting or other markings on the paper.

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below _____

HOLDING TIMES

The objective of this parameter is to ascertain the validity of the results based on the holding time of the sample from time of collection to the time of analysis.

Complete table for all samples and note the analysis and/or preservation not within criteria

SAMPLE ID	DATE SAMPLED	DATE EXTRACTED/ANALYZED	pH	ACTION
All samples extracted and analyzed within method recommended holding time. Sample preservation was acceptable.				

Cooler temperature (Criteria: $4 \pm 2^{\circ}\text{C}$): 3.6°C

Actions

Results will be qualified based on the criteria of the following Table:

Table 1. Holding Time Actions for Semivolatile Analyses

Matrix	Preserved	Criteria	Action	
			Detected Associated Compounds	Non-Detected Associated Compounds
Aqueous	No	≤ 7 days (for extraction) ≤ 40 days (for analysis)	Use professional judgment	
	No	> 7 days (for extraction) > 40 days (for analysis)	J	Use professional judgment
	Yes	≤ 7 days (for extraction) ≤ 40 days (for analysis)	No qualification	
	Yes	> 7 days (for extraction) > 40 days (for analysis)	J	UJ
	Yes/No	Grossly Exceeded	J	UJ or R
Non-Aqueous	No	≤ 14 days (for extraction) ≤ 40 days (for analysis)	Use professional judgment	
	No	> 14 days (for extraction) > 40 days (for analysis)	J	Use professional judgment
	Yes	≤ 14 days (for extraction) ≤ 40 days (for analysis)	No qualification	
	Yes	> 14 days (for extraction) > 40 days (for analysis)	J	UJ
	Yes/No	Grossly Exceeded	J	UJ or R

DATA REVIEW WORKSHEETS

All criteria were met X
Criteria were not met see below _____

GC/MS TUNING

The assessment of the tuning results is to determine if the sample instrumentation is within the standard tuning QC limits

X The DFTPP performance results were reviewed and found to be within the specified criteria.

X DFTPP tuning was performed for every 12 hours of sample analysis.

If no, use professional judgment to determine whether the associated data should be accepted, qualified or rejected.

Notes: These requirements do not apply when samples are analyzed by the Selected Ion Monitoring (SIM) technique.

All mass spectrometer conditions must be identical to those used during the sample analysis. Background subtraction actions resulting in spectral distortion are unacceptable

Notes: No data should be qualified based of DFTPP failure.

The requirement to analyze the instrument performance check solution is optional when analysis of PAHs/pentachlorophenol is to be performed by the SIM technique.

List the samples affected:

[illegible]

Actions:

1. If sample are analyzed without a preceding valid instrument performance check or are analyzed 12 hours after the Instrument Performance Check, qualify all data in those samples as unusable (R).
2. If ion abundance criteria are not met, use professional judgment to determine to what extent the data may be utilized.
3. State in the Data Review Narrative, decisions to use analytical data associated with DFTPP instrument performance checks not meeting the contract requirements.
4. Use professional judgment to determine if associated data should be qualified based on the spectrum of the mass calibration compounds.

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

INITIAL CALIBRATION VERIFICATION

Compliance requirements for satisfactory instrument calibration are established to ensure that the instrument is capable of producing and maintaining acceptable quantitative data.

Date of initial calibration: 04/21/2016 (SIM) 04/04-05/2016 (Scan)
 Instrument ID numbers: GCMS3M GCMSF
 Matrix/Level: Aqueous/low Aqueous/low
 Date of initial calibration: 04/26/16; 04/27/16 (Scan) 04/13-14/16 (Scan)
 Instrument ID numbers: GCMS3E GCMSZ
 Matrix/Level: Aqueous/low

DATE	LAB FILE ID#	CRITERIA OUT RFs, %RSD, %D, r	COMPOUND	SAMPLES AFFECTED
Initial and initial calibration verification meets the method and guidance validation document performance criteria.				

Actions:

Qualify the initial calibration analytes listed in Table 2 using the following criteria:

Table 3. Initial Calibration Actions for Semivolatile Analysis

Criteria	Action	
	Detect	Non-detect
Initial Calibration not performed at specified frequency and sequence	Use professional judgment R	Use professional judgment R
Initial Calibration not performed at the specified concentrations	J	UJ
RRF < Minimum RRF in Table 2 for target analyte	Use professional judgment J+ or R	R
RRF ≥ Minimum RRF in Table 2 for target analyte	No qualification	No qualification
%RSD > Maximum %RSD in Table 2 for target analyte	J	Use professional judgment
%RSD ≤ Maximum %RSD in Table 2 for target analyte	No qualification	No qualification

Initial Calibration**Table 2. RRF, %RSD, and %D Acceptance Criteria in Initial Calibration and CCV for Semivolatile Analysis**

Analyte	Minimum RRF	Maximum %RSD	Opening Maximum %D ¹	Opening Maximum %D ¹
1,4-Dioxane	0.010	40.0	± 40.0	± 50.0
Benzaldehyde	0.100	40.0	± 40.0	± 50.0
Phenol	0.080	20.0	± 20.0	± 25.0
Bis(2-chloroethyl)ether	0.100	20.0	± 20.0	± 25.0
2-Chlorophenol	0.200	20.0	± 20.0	± 25.0
2-Methylphenol	0.010	20.0	± 20.0	± 25.0
3-Methylphenol	0.010	20.0	± 20.0	± 25.0
2,2'-Oxybis-(1-chloropropane)	0.010	20.0	± 25.0	± 50.0
Acetophenone	0.060	20.0	± 20.0	± 25.0
4-Methylphenol	0.010	20.0	± 20.0	± 25.0
N-Nitroso-di-n-propylamine	0.080	20.0	± 25.0	± 25.0
Hexachloroethane	0.100	20.0	± 20.0	± 25.0
Nitrobenzene	0.090	20.0	± 20.0	± 25.0
Isophorone	0.100	20.0	± 20.0	± 25.0
2-Nitrophenol	0.060	20.0	± 20.0	± 25.0
2,4-Dimethylphenol	0.050	20.0	± 25.0	± 50.0
Bis(2-chloroethoxy)methane	0.080	20.0	± 20.0	± 25.0
2,4-Dichlorophenol	0.060	20.0	± 20.0	± 25.0
Naphthalene	0.200	20.0	± 20.0	± 25.0
4-Chloroaniline	0.010	40.0	± 40.0	± 50.0
Hexachlorobutadiene	0.040	20.0	± 20.0	± 25.0
Caprolactam	0.010	40.0	± 30.0	± 50.0
4-Chloro-3-methylphenol	0.040	20.0	± 20.0	± 25.0
2-Methylnaphthalene	0.100	20.0	± 20.0	± 25.0
Hexachlorocyclopentadiene	0.010	40.0	± 40.0	± 50.0
2,4,6-Trichlorophenol	0.090	20.0	± 20.0	± 25.0
2,4,5-Trichlorophenol	0.100	20.0	± 20.0	± 25.0
1,1'-Biphenyl	0.200	20.0	± 20.0	± 25.0

DATA REVIEW WORKSHEETS

Analyte	Minimum RRF	Maximum %RSD	Opening Maximum %D ¹	Opening Maximum %D ¹
2-Chloronaphthalene	0.300	20.0	± 20.0	± 25.0
2-Nitroaniline	0.060	20.0	± 25.0	± 25.0
Dimethylphthalate	0.300	20.0	± 25.0	± 25.0
2,6-Dinitrotoluene	0.080	20.0	± 20.0	± 25.0
Acenaphthylene	0.400	20.0	± 20.0	± 25.0
3-Nitroaniline	0.010	20.0	± 25.0	± 50.0
Acenaphthene	0.200	20.0	± 20.0	± 25.0
2,4-Dinitrophenol	0.010	40.0	± 50.0	± 50.0
4-Nitrophenol	0.010	40.0	± 40.0	± 50.0
Dibenzofuran	0.300	20.0	± 20.0	± 25.0
2,4-Dinitrotoluene	0.070	20.0	± 20.0	± 25.0
Diethylphthalate	0.300	20.0	± 20.0	± 25.0
1,2,4,5-Tetrachlorobenzene	0.100	20.0	± 20.0	± 25.0
4-Chlorophenyl-phenylether	0.100	20.0	± 20.0	± 25.0
Fluorene	0.200	20.0	± 20.0	± 25.0
4-Nitroaniline	0.010	40.0	± 40.0	± 50.0
4,6-Dinitro-2-methylphenol	0.010	40.0	± 30.0	± 50.0
4-Bromophenyl-phenyl ether	0.070	20.0	± 20.0	± 25.0
N-Nitrosodiphenylamine	0.100	20.0	± 20.0	± 25.0
Hexachlorobenzene	0.050	20.0	± 20.0	± 25.0
Atrazine	0.010	40.0	± 25.0	± 50.0
Pentachlorophenol	0.010	40.0	± 40.0	± 50.0
Phenanthrene	0.200	20.0	± 20.0	± 25.0
Anthracene	0.200	20.0	± 20.0	± 25.0
Carbazole	0.050	20.0	± 20.0	± 25.0
Di-n-butylphthalate	0.500	20.0	± 20.0	± 25.0
Fluoranthene	0.100	20.0	± 20.0	± 25.0
Pyrene	0.400	20.0	± 25.0	± 50.0
Butylbenzylphthalate	0.100	20.0	± 25.0	± 50.0

DATA REVIEW WORKSHEETS

Analyte	Minimum RRF	Maximum %RSD	Opening Maximum %D ¹	Opening Maximum %D ¹
3,3'-Dichlorobenzidine	0.010	40.0	± 40.0	± 50.0
Benzo(a)anthracene	0.300	20.0	± 20.0	± 25.0
Chrysene	0.200	20.0	± 20.0	± 50.0
Bis(2-ethylhexyl) phthalate	0.200	20.0	± 25.0	± 50.0
Di-n-octylphthalate	0.010	40.0	± 40.0	± 50.0
Benzo(b)fluoranthene	0.010	20.0	± 25.0	± 50.0
Benzo(k)fluoranthene	0.010	20.0	± 25.0	± 50.0
Benzo(a)pyrene	0.010	20.0	± 20.0	± 50.0
Indeno(1,2,3-cd)pyrene	0.010	20.0	± 25.0	± 50.0
Dibenzo(a,h)anthracene	0.010	20.0	± 25.0	± 50.0
Benzo(g,h,i)perylene	0.010	20.0	± 30.0	± 50.0
2,3,4,6-Tetrachlorophenol	0.040	20.0	± 20.0	± 50.0
Naphthalene	0.600	20.0	± 25.0	± 25.0
2-Methylnaphthalene	0.300	20.0	± 20.0	± 25.0
Acenaphthylene	0.900	20.0	± 20.0	± 25.0
Acenaphthene	0.500	20.0	± 20.0	± 25.0
Fluorene	0.700	20.0	± 25.0	± 50.0
Phenanthrene	0.300	20.0	± 25.0	± 50.0
Anthracene	0.400	20.0	± 25.0	± 50.0
Fluoranthene	0.400	20.0	± 25.0	± 50.0
Pyrene	0.500	20.0	± 30.0	± 50.0
Benzo(a)anthracene	0.400	20.0	± 25.0	± 50.0
Chrysene	0.400	20.0	± 25.0	± 50.0
Benzo(b)fluoranthene	0.100	20.0	± 30.0	± 50.0
Benzo(k)fluoranthene	0.100	20.0	± 30.0	± 50.0
Benzo(a)pyrene	0.100	20.0	± 25.0	± 50.0
Indeno(1,2,3-cd)pyrene	0.100	20.0	± 40.0	± 50.0
Dibenzo(a,h)anthracene	0.010	25.0	± 40.0	± 50.0
Benzo(g,h,i)perylene	0.020	25.0	± 40.0	± 50.0

DATA REVIEW WORKSHEETS

Pentachlorophenol	0.010	40.0	- 50.0	+ 50.0
Deuterated Monitoring Compounds				
Analyte	Minimum RRF	Maximum %RSD	Opening Maximum %D ¹	Closing Maximum %D
1,4-Dioxane-d ₈	0.010	20.0	- 25.0	± 50.0
Phenol-d ₅	0.010	20.0	- 25.0	± 25.0
Bis-(2-chloroethyl)ether-d ₈	0.100	20.0	- 20.0	+ 25.0
2-Chlorophenol-d ₁	0.200	20.0	- 20.0	± 25.0
4-Methylphenol-d ₈	0.010	20.0	- 20.0	+ 25.0
4-Chloroaniline-d ₄	0.010	40.0	± 40.0	± 50.0
Nitrobenzene-d ₅	0.050	20.0	± 20.0	± 25.0
2-Nitrophenol-d ₄	0.050	20.0	- 20.0	± 25.0
2,4-Dichlorophenol-d ₃	0.060	20.0	- 20.0	+ 25.0
Dimethylphthalate-d ₆	0.300	20.0	- 20.0	+ 25.0
Acenaphthylene-d ₈	0.400	20.0	- 20.0	+ 25.0
4-Nitrophenol-d ₄	0.010	40.0	- 40.0	± 50.0
Fluorene-d ₁₀	0.100	20.0	± 20.0	± 25.0
4,6-Dinitro-2-methylphenol-d ₇	0.010	40.0	- 30.0	± 50.0
Anthracene-d ₁₀	0.300	20.0	- 20.0	+ 25.0
Pyrene-d ₁₀	0.300	20.0	- 25.0	± 50.0
Benzo(a)pyrene-d ₁₂	0.010	20.0	- 20.0	+ 50.0
Fluoranthene-d ₁₀ (SIM)	0.400	20.0	- 25.0	± 50.0
2-Methylnaphthalene-d ₁₀ (SIM)	0.300	20.0	- 20.0	+ 25.0

¹ If a closing CCV is acting as an opening CCV, all target analytes must meet the requirements for an opening CCV.

Note: If analysis by SIM technique is requested for PAH/pentachlorophenols, calibration standards analyzed at 0.10, 0.20, 0.40, 0.80, and 1.0 ng/uL for each target compound of interest and the associated DMCs. Pentachlorophenol will require only a four point initial calibration at 0.20, 0.40, 0.80, and 1.0 ng/uL.

DATA REVIEW WORKSHEETS

All criteria were met _____
 Criteria were not met _____
 and/or see below X

CONTINUING CALIBRATION VERIFICATION

Compliance requirements for satisfactory instrument calibration are established to ensure that the instrument is capable of producing and maintaining acceptable quantitative data.

Date of initial calibration: 04/04-05/16 (Scan)
 Date of initial calibration verification (ICV): 04/05-06/16
 Date of continuing calibration verification (CCV): 05/03/16
 Date of closing CCV: _____
 Instrument ID numbers: GCMSF
 Matrix/Level: Aqueous/low

Date of initial calibration: 04/21/16 (SIM)
 Date of initial calibration verification (ICV): 04/21/16
 Date of continuing calibration verification (CCV): 05/04-05/16; 05/06/16
 Date of closing CCV: _____
 Instrument ID numbers: GCMS3M
 Matrix/Level: Aqueous/low

Date of initial calibration: 04/26/16; 04/27/16 (Scan)
 Date of initial calibration verification (ICV): 04/26-27/16; 04/27/16
 Date of continuing calibration verification (CCV): 05/11/16
 Date of closing CCV: _____
 Instrument ID numbers: GCMS3E
 Matrix/Level: Aqueous/low

Date of initial calibration: 04/13-14/16 (Scan)
 Date of initial calibration verification (ICV): 04/14/16
 Date of continuing calibration verification (CCV): 05/05/16
 Date of closing CCV: _____
 Instrument ID numbers: GCMS3E
 Matrix/Level: Aqueous/low

DATE	LAB FILE ID#	CRITERIA OUT RFs, %RSD, %D, r	COMPOUND	SAMPLES AFFECTED
GCMS3E				
05/11/16	cc3621-21	-22.2	1,2,4,5-tetrachlorobenzene	JC19473-3
		-24.5	Atrazine*	
GCMSF				
05/03/16	cc6563-25	-24.4/-23.4	n-Nitroso-di-n-propylamine	JC19473-1; -2
		-20.7/-25.9	2-Nitroaniline	
		29.3	4-Nitrophenol*	
05/03/16	cc6564-25	-20.8/-22.4	1,2,4,5-tetrachlorobenzene	JC19473-1; -2

DATA REVIEW WORKSHEETS

Note: Initial and continuing calibration verifications meet the method and guidance document required performance criteria except the cases described above. Results qualified as estimated (J), (UJ) for non-detects.

No closing calibration verification included in data package. No action taken, professional judgment.

* Analytes with % difference in the continue calibration verification outside the method performance criteria but within the validation guidelines criteria, + 40 %. No action taken.

QC samples for this data package were analyzed in the instrument GCMSZ (Scan). Analytes not meeting continue calibration verification criteria were not qualified.

Actions:

Notes: Verify that the CCV is run at the required frequency (an opening and closing CCV must be run within 12-hour period).

All DMCs must meet the RRF values given in Table 2. No qualification of the data is necessary on DMCs RRF and %RSD/%D alone. Use professional judgment to evaluate DMCs and %RSD/%D data in conjunction with DMCs recoveries to determine the need for qualification of the data.

Qualify the initial calibration analytes listed in Table 2 using the following criteria in the CCVs:

Table 4. CCV Actions for Semivolatile Analysis

Criteria for Opening CCV	Criteria for Closing CCV	Action	
		Detect	Non-detect
CCV not performed at required frequency and sequence	CCV not performed at required frequency	Use professional judgment R	Use professional judgment R
CCV not performed at specified concentration	CCV not performed at specified concentration	Use professional judgment	Use professional judgment
RRF < Minimum RRF in Table 2 for target analyte	RRF < Minimum RRF in Table 2 for target analyte	Use professional judgment J or R	R
RRF ≥ Minimum RRF in Table 2 for target analyte	RRF ≥ Minimum RRF in Table 2 for target analyte	No qualification	No qualification
%D outside the Opening Maximum %D limits in Table 2 for target analyte	%D outside the Closing Maximum %D limits in Table 2 for target analyte	J	UJ
%D within the inclusive Opening Maximum %D limits in Table 2 for target analyte	%D within the inclusive Closing Maximum %D limits in Table 2 for target analyte	No qualification	No qualification

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

BLANK ANALYSIS RESULTS (Sections 1 & 2)

The assessment of the blank analysis results is to determine the existence and magnitude of contamination problems. The criteria for evaluation of blanks apply only to blanks associated with the samples, including trip, equipment, and laboratory blanks. If problems with any blanks exist, all data associated with the case must be carefully evaluated to determine whether or not there is an inherent variability in the data for the case, or if the problem is an isolated occurrence not affecting other data.

List the contamination in the blanks below. High and low levels blanks must be treated separately.

Notes: The concentration of non-target compounds in all blanks must be less than or equal to 10 ug/L.

The concentration of target compounds in all blanks must be less than its CRQL listed in the method.

Samples taken from a drinking water tap do not have and associated field blank.

Laboratory blanks

DATE ANALYZED	LAB ID	LEVEL/ MATRIX	COMPOUND	CONCENTRATION UNITS
_No_target_analytes_detected_in_method_blanks._				

Field/Equipment/Trip blank

DATE ANALYZED	LAB ID	LEVEL/ MATRIX	COMPOUND	CONCENTRATION UNITS
_No_target_analytes_detected_in_the_equipment_blank._No_field/trip_blanks_analyzed_with_this_data_package._				

All criteria were met X
 Criteria were not met
 and/or see below

BLANK ANALYSIS RESULTS (Section 3)

Blank Actions

Qualify samples based on the criteria summarized in Table 5:

Table 5. Blank and TCLP/SPLP LEB Actions for Semivolatile Analysis

Blank Type	Blank Result	Sample Result	Action
Method, TCLP/SPLP LEB, Field	Detect	Non-detect	No qualification
	< CRQL	< CRQL	Report at CRQL and qualify as non-detect (U)
		≥ CRQL	Use professional judgment
	≥ CRQL	< CRQL	Report at CRQL and qualify as non-detect (U)
		≥ CRQL but < Blank Result	Report at sample results and qualify as non-detect (U) or as unusable (R)
		≥ CRQL and ≥ Blank Result	Use professional judgment
	Grossly high	Detect	Report at sample results and qualify as unusable (R)
	TIC > 5.0 ug/L (water) or 0.0050 mg/L (TCLP leachate) or TIC > 170 ug/Kg (soil)	Detect	Use professional judgment

List samples qualified

CONTAMINATION SOURCE/LEVEL	COMPOUND	CONC/UNITS	AL/UNITS	SQL	AFFECTED SAMPLES

All criteria were met X
 Criteria were not met
 and/or see below

SURROGATE SPIKE RECOVERIES – DEUTERATED MONITORING COMPOUNDS (DMCs)

Laboratory performance of individual samples is established by evaluation of surrogate spike recoveries – deuterated monitoring compounds. All samples are spiked with surrogate compounds prior to sample analysis. The accuracy of the analysis is measured by the surrogate percent recovery. Since the effects of the sample matrix are frequently outside the control of the laboratory and may present relatively unique problems, the validation of data is frequently subjective and demands analytical experience and professional judgment.

Notes: Recoveries for DMCs in samples and blanks must be within the limits specified in Table 6.

The recovery limits for any of the compounds listed in Table 6 may be expanded at any time during the period of performance if USEPA determines that the limits are too restrictive.

If a DMC is not added in the samples and blanks or the concentrations of DMCs in the samples and blank not the specified, use professional judgment in qualifying the data.

Table 7. DMC Actions for Semivolatile Analysis

Criteria	Action	
	Detect	Non-detect
%R < 10% (excluding DMCs with 10% as a lower acceptance limit)	J-	R
10% ≤ %R (excluding DMCs with 10% as a lower acceptance limit) < Lower Acceptance Limit	J-	UJ
Lower Acceptance limit ≤ %R ≤ Upper Acceptance Limit	No qualification	No qualification
%R > Upper Acceptance Limit	J+	No qualification

List the percent recoveries (%Rs) which do not meet the criteria for DMCs (surrogate) recovery.

Matrix: Groundwater

SAMPLE ID	SURROGATE COMPOUND	ACTION
-----------	--------------------	--------

DMCs meet the required criteria. Non-deuterated surrogates added to the samples were
within laboratory recovery limits.

Table 8. Semivolatile DMCs and the Associated Target Analytes

1,4-Dioxane-d₈ (DMC-1)	Phenol-d₅ (DMC-2)	Bis(2-Chloroethyl) ether-d₈ (DMC-3)
1,4-Dioxane	Benzaldehyde Phenol	Bis(2-chloroethyl) ether 2,2'-Oxybis(1-chloropropane) Bis(2-chloroethoxy) methane
2-Chlorophenol-d₄ (DMC-4)	4-Methylphenol-d₄ (DMC-5)	4-Chloroaniline-d₄ (DMC-6)
2-Chlorophenol	2-Methylphenol 3-Methylphenol 4-Methylphenol 2,4-Dimethylphenol	4-Chloroaniline Hexachlorocyclopentadiene Dichlorobenzidine
Nitrobenzene-d₅ (DMC-7)	2-Nitrophenol-d₄ (DMC-8)	2,4-Dichlorophenol-d₃ (DMC-9)
Acetophenone N-Nitroso-di-n-propylamine Hexachloroethane Nitrobenzene 2,6-Dinitrotoluene 2,4-Dinitrotoluene N-Nitrosodiphenylamine	Isophorone 2-Nitrophenol	2,4-Dichlorophenol Hexachlorobutadiene Hexachlorocyclopentadiene 4-Chloro-3-methylphenol 2,4,6-Trichlorophenol 2,4,5-Trichlorophenol 1,2,4,5-Tetrachlorobenzene *Pentachlorophenol 2,3,4,6-Tetrachlorophenol
Dimethylphthalate-d₈ (DMC-10)	Acenaphthylene-d₈ (DMC-11)	4-Nitrophenol-d₄ (DMC-12)
Caprolactam 1,1'-Biphenyl Dimethylphthalate Diethylphthalate Di-n-butylphthalate Butylbenzylphthalate Bis(2-ethylhexyl) phthalate Di-n-octylphthalate	*Naphthalene *2-Methylnaphthalene 2-Chloronaphthalene *Acenaphthylene *Acenaphthene	2-Nitroaniline 3-Nitroaniline 2,4-Dinitrophenol 4-Nitrophenol 4-Nitroaniline

DATA REVIEW WORKSHEETS

Fluorene-d₁₀ (DMC-13)	4,6-Dinitro-2-methylphenol-d₂ (DMC-14)	Anthracene-d₁₀ (DMC-15)
Dibenzofuran *Fluorene 4-Chlorophenyl-phenylether 4-Bromophenyl-phenylether Carbazole	4,6-Dinitro-2-methylphenol	Hexachlorobenzene Atrazine *Phenanthrene *Anthracene
Pyrene-d₁₀ (DMC-16)	Benzo(a)pyrene-d₁₂ (DMC-17)	
*Fluoranthene *Pyrene *Benzo(a)anthracene *Chrysene	3,3'-Dichlorobenzidine *Benzo(b)fluoranthene *Benzo(k)fluoranthene *Benzo(a)pyrene *Indeno(1,2,3-cd)pyrene *Dibenzo(a,h)anthracene *Benzo(g,h,i)perylene	

*Included in optional Target Analyte List (TAL) of PAHs and PCP only.

Table 9. Semivolatile SIM DMCs and the Associated Target Analytes

Fluoranthene-d₁₀ (DMC-1)	2-Methylnaphthalene-d₁₀ (DMC-2)
Fluoranthene	Naphthalene
Pyrene	2-Methylnaphthalene
Benzo(a)anthracene	Acenaphthylene
Chrysene	Acenaphthene
Benzo(b)fluoranthene	Fluorene
Benzo(k)fluoranthene	Pentachlorophenol
Benzo(a)pyrene	Phenanthrene
Indeno(1,2,3-cd)pyrene	Anthracene
Dibenzo(a,h)anthracene	
Benzo(g,h,i)perylene	

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met _____
 and/or see below _____

VII. A MATRIX SPIKE/MATRIX SPIKE DUPLICATE (MS/MSD)

This data is generated to determine long term precision and accuracy in the analytical method for various matrices. This data alone cannot be used to evaluate the precision and accuracy of individual samples. If any % R in the MS or MSD falls outside the designated range, the reviewer should determine if there are matrix effects, i.e. LCS data are within the QC limits but MS/MSD data are outside QC limit.

1. MS/MSD Recoveries and Precision Criteria

The laboratory should use one MS and a duplicate analysis of an unspiked field sample if target analytes are expected in the sample. If target analytes are not expected, MS/MSD should be analyzed.

NOTES: Data for MS and MSDs will not be present unless requested by the Region.
 Notify the Contract Laboratory COR if a field or trip blank was used for the MS and MSD.

For a Matrix Spike that does not meet criteria, apply the action to only the field sample used to prepare the Matrix Spike sample. If it is clearly stated in the data validation materials that the samples were taken through incremental sampling or some other method guaranteeing the homogeneity of the sample group, then the entire sample group may be qualified.

List the %Rs, RPD of the compounds which do not meet the criteria.

Sample ID: <u> JC19295-1* </u>	Matrix/Level: <u> Groundwater </u>
Sample ID: <u> JC19332-3_(SIM) </u>	Matrix/Level: <u> Groundwater </u>
Sample ID: <u> JC19414-1 </u>	Matrix/Level: <u> Soil </u>
Sample ID: <u> JC19423-1_(SIM) </u>	Matrix/Level: <u> Soil </u>

* see enclosed list

MS OR MSD	COMPOUND	% R	RPD	QC LIMITS	ACTION
<u> JC19423-1MS/MSD </u>					
<u> MS/MSD </u>	<u> 1,4-Dioxane </u>	<u> 10/10 </u> %		<u> 50 - 150 </u>	<u> No action </u>

Note: No action taken, MS/MSD results apply only to unspiked sample.

- * QC limits are laboratory in-house performance criteria, LL = lower limit, UL = upper limit.
- * If QC limits are not available, use limits of 70 – 130 %.

DATA REVIEW WORKSHEET

MATRIX SPIKE/MATRIX SPIKE DUPLICATE

The QC reported here applies to the following samples:
JC19473-1, JC19473-2

Method: SW846 8270D

Compound	JC19295-1 ug/l	Q	Spike ug/l	MS ug/l	MS %	Spike ug/l	MSD ug/l	MSD %	RPD	Limits Rec/RPD
2-Chlorophenol	ND		100	74.1	74	100	95.8	96	26* a	49-110/20
4-Chloro-3-methyl phenol	ND		100	74.3	74	100	89.5	90	19* a	44-121/18
2-Methylphenol	ND		100	82.9	83	100	104	104	23* a	47-112/18
3&4-Methylphenol	ND		100	82.1	82	100	103	103	23* a	44-113/19
Phenol	ND		100	64.2	64	100	83.1	83	26* a	22-100/22
2,4,5-Trichlorophenol	ND		100	74.0	74	100	92.1	92	22* a	51-124/20
2,4,6-Trichlorophenol	ND		100	79.7	80	100	99.8	100	22* a	53-120/21
Acetophenone	ND		100	236	236* b	100	279	279* b	17	31-141/23
Benzaldehyde	ND		100	385	385* b	100	435	435* b	12	10-164/30
1-Methylnaphthalene	ND		100	126	126* b	100	141	141* b	11	34-124/25

Note: No action taken, MS/MSD results apply to unspiked sample.

DATA REVIEW WORKSHEETS

Actions:

QUALITY	%R < LL	%R > UL
Positive results	J	J
Nondetects results	R	Accept

MS/MSD criteria apply only to the unspiked sample, its dilutions, and the associated MS/MSD samples:

If the % R for the affected compounds were < LL (or 70 %), qualify positive results (J) and nondetects (JJ).

If the % R for the affected compounds were > UL (or 130 %), only qualify positive results (J).

If 25 % or more of all MS/MSD %R were < LL (or 70 %) or if two or more MS/MSD %Rs were < 10%, qualify all positive results (J) and reject nondetects (R).

A separate worksheet should be used for each MS/MSD pair.

All criteria were met ☒ X
Criteria were not met
and/or see below _____

INTERNAL STANDARD PERFORMANCE

The assessment of the internal standard (IS) parameter is used to assist the data reviewer in determining the condition of the analytical instrumentation.

List the internal standard area of samples which do not meet the criteria.

DATE	SAMPLE ID	IS OUT	IS AREA	ACCEPTABLE RANGE	ACTION
------	-----------	--------	---------	------------------	--------

Internal standard area counts meet the required criteria.

Action:

1. If an internal standard area count for a sample or blank is greater than 200.0% of the area for the associated standard (opening CCV or mid-point standard from initial calibration) (see Table 10 below):
 - a. Qualify detects for compounds quantitated using that internal standard as estimated low (J-).
 - b. Do not qualify non-detected associated compounds.
2. If an internal standard area count for a sample or blank is less than 20.0% of the area for the associated standard (opening CCV or mid-point standard from initial calibration):
 - a. Qualify detects for compounds quantitated using that internal standard as estimated high (J+).
 - b. Qualify non-detected associated compounds as unusable (R).
3. If an internal standard area count for a sample or blank is greater than or equal to 50.0%, and less than or equal to 200% of the area for the associated standard opening CCV or mid-point standard from initial calibration, no qualification of the data is necessary.
4. If an internal standard RT varies by more than 10.0 seconds: Examine the chromatographic profile for that sample to determine if any false positives or negatives exist. For shifts of a large magnitude, the reviewer may consider partial or total rejection of the data for that sample fraction. Detects should not need to be qualified as unusable (R) if the mass spectral criteria are met.
5. If an internal standard RT varies by less than or equal to 10.0 seconds, no qualification of the data is necessary.

Note: Inform the Contract Laboratory Program Project Officer (CLP PO) if the internal standard performance criteria are grossly exceeded. Note in the Data Review Narrative potential effects on the data resulting from unacceptable internal standard performance.

DATA REVIEW WORKSHEETS

State in the Data Review Narrative if the required internal standard compounds are not added to a sample or blank or if the required internal standard compound is not analyzed at the specified concentration.

Actions:

Table 10. Internal Standard Actions for Semivolatile Analysis

Criteria	Action	
	Detect	Non-detect
Area response < 20% of the opening CCV or mid-point standard CS3 from ICAL	J+	R
20% ≤ Area response < 50% of the opening CCV or mid-point standard CS3 from ICAL	J+	UJ
50% ≤ Area response ≤ 200% of the opening CCV or mid-point standard CS3 from ICAL	No qualification	No qualification
Area response > 200% of the opening CCV or mid-point standard CS3 from ICAL	J-	No qualification
RT shift between sample/blank and opening CCV or mid-point standard CS3 from ICAL > 10.0 seconds	R	R
RT shift between sample/blank and opening CCV or mid-point standard CS3 from ICAL < 10.0 seconds	No qualification	No qualification

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

TARGET COMPOUND IDENTIFICATION

Criteria:

Is the Relative Retention Times (RRTs) of reported compounds within ± 0.06 RRT units of the standard RRT [opening Continuing Calibration Verification (CCV) or mid-point standard from the initial calibration]. Yes? or No?

List compounds not meeting the criteria described above:

Sample ID	Compounds	Actions
=====	=====	=====
_____	_____	_____
_____	_____	_____
_____	_____	_____

Mass spectra of the sample compound and a current laboratory-generated standard [i.e., the mass spectrum from the associated calibration standard (opening CCV or mid-point standard from initial calibration)] must match according to the following criteria:

- a. All ions present in the standard mass spectrum at a relative intensity greater than 10% must be present in the sample spectrum.
- b. The relative intensities of these ions must agree within $\pm 20\%$ between the standard and sample spectra (e.g., for an ion with an abundance of 50% in the standard spectrum, the corresponding sample ion abundance must be between 30-70%).
- c. Ions present at greater than 10% in the sample mass spectrum, but not present in the standard spectrum, must be evaluated by a reviewer experienced in mass spectral interpretation.

List compounds not meeting the criteria described above:

Sample ID	Compounds	Actions
=====	=====	=====
_____	_____	_____
_____	_____	_____
_____	_____	_____

Identified compounds meet the required criteria

DATA REVIEW WORKSHEETS

Action:

1. The application of qualitative criteria for GC/MS analysis of target compounds requires professional judgment. It is up to the reviewer's discretion to obtain additional information from the laboratory. If it is determined that incorrect identifications were made, qualify all such data as unusable (R).
2. Use professional judgment to qualify the data if it is determined that cross-contamination has occurred.
3. Note in the Data Review Narrative any changes made to the reported compounds or concerns regarding target compound identifications. Note, for Contract Laboratory COR action, the necessity for numerous or significant changes.

TENTATIVELY IDENTIFIED COMPOUNDS (TICS)

NOTE: Tentatively identified compounds should only be evaluated when requested by a party from outside of the Hazardous Waste Support Section (HWSS).

List TICs

Sample ID	Compound	Sample ID	Compound
=====			

Action:

1. Qualify all TIC results for which there is presumptive evidence of a match (e.g. greater than or equal to 85% match) as tentatively identified (NJ), with approximated concentrations. TICs labeled "unknown" are qualified as estimated (J).
2. General actions related to the review of TIC results are as follows:
 - a. If it is determined that a tentative identification of a non-target compound is unacceptable, change the tentative identification to "unknown" or another appropriate identification, and qualify the result as estimated (J).
 - b. If all contractually-required peaks were not library searched and quantitated, the Region's designated representative may request these data from the laboratory.
3. In deciding whether a library search result for a TIC represents a reasonable identification, use professional judgment. If there is more than one possible match, report the result as "either compound X or compound Y". If there is a lack of isomer specificity, change the TIC result to a nonspecific isomer result (e.g., 1,3,5-trimethyl benzene to trimethyl benzene isomer) or to a compound class (e.g., 2-methyl, 3-ethyl benzene to a substituted aromatic compound).
4. The reviewer may elect to report all similar compounds as a total (e.g., all alkanes may be summarized and reported as total hydrocarbons).

DATA REVIEW WORKSHEETS

5. Target compounds from other fractions and suspected laboratory contaminants should be marked as "non-reportable".
6. Other Case factors may influence TIC judgments. If a sample TIC match is poor, but other samples have a TIC with a valid library match, similar RRT, and the same ions, infer identification information from the other sample TIC results.
7. Note in the Data Review Narrative any changes made to the reported data or any concerns regarding TIC identifications.
8. Note, for Contract Laboratory COR action, failure to properly evaluate and report TICs

All criteria were met X
 Criteria were not met
 and/or see below

SAMPLE QUANTITATION AND REPORTED CONTRACT REQUIRED QUANTITATION LIMITS (CRQLS)

Action:

1. When a sample is analyzed at more than one dilution, the lower CRQL are used unless a QC exceedance dictates the use of higher CRQLs from the diluted sample. Samples reported with an "E" qualifier should be reported from the diluted sample.
2. If any discrepancies are found, the Region's designated representative may contact the laboratory to obtain additional information that could resolve any differences. If a discrepancy remains unresolved, the reviewer must use professional judgment to decide which value is the most accurate. Under these circumstances, the reviewer may determine that qualification of data is warranted. Note in the Data Review Narrative a description of the reasons for data qualification and the qualification that is applied to the data.
3. For non-aqueous samples, if the solids is less than 10.0%, use professional judgment for both detects and non-detects. If the percent solid for a soil sample is greater than or equal to 10.0% and less than 30.0%, use professional judgment to qualify detects and non-detects. If the percent solid for a soil sample is greater than or equal to 30.0%, detects and non-detects should not be qualified (see Table 11).
4. Note, for Contract Laboratory COR action, numerous or significant failures to accurately quantify the target compounds or to properly evaluate and adjust CRQLs.
5. Results between MDL and CRQL should be qualified as estimated "J".
6. Results < MDL should be reported at the CRQL and qualified "U". MDLs themselves should not be reported.

Table 11. Percent Solids Actions for Semivolatile Analysis for Non-Aqueous Samples

Criteria	Action	
	Detects	Non-detects
%Solids < 10.0%	Use professional judgment	Use professional judgment
10.0% ≤ %Solids ≤ 30.0%	Use professional judgment	Use professional judgment
%Solids > 30.0%	No qualification	No qualification

SAMPLE QUANTITATION

The sample quantitation evaluation is to verify laboratory quantitation results. In the space below, please show a minimum of one sample calculation:

Sample ID: JC19473-2 Analyte: 1,4-Dioxane RF: 0.207

$$\begin{aligned}
 [] &= (41714)(4.0)/(341879)(0.207) \\
 &= 2.36 \text{ ppm} \quad \text{Ok}
 \end{aligned}$$

DATA REVIEW WORKSHEETS

QUANTITATION LIMITS

A. Dilution performed

SAMPLE ID	DILUTION FACTOR	REASON FOR DILUTION

DATA REVIEW WORKSHEETS

All criteria were met N/A
 Criteria were not met
 and/or see below

FIELD DUPLICATE PRECISION

Sample IDs: -

Matrix: -

Field duplicate samples may be taken and analyzed as an indication of overall precision. These analyses measure both field and lab precision; therefore, the results may have more variability than laboratory duplicates which only laboratory performance. It is also expected that soil duplicate results will have a greater variance than water matrices due to difficulties associated with collecting identical field duplicate samples.

The project QAPP should be reviewed for project-specific information.

Suggested criteria: if large RPD (> 50 %) is observed, confirm identification of the samples and note differences. If both samples and duplicate are <5 SQL, the RPD criteria is doubled.

COMPOUND	SQL ug/L	SAMPLE CONC.	DUPLICATE CONC.	RPD	ACTION
No field/laboratory duplicate analyzed as part of this data package. MS/MSD % and blank spike/blank spike duplicate recoveries RPD used to assess precision; RPD within the required criteria < 50 % for detected target analytes.					

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

OTHER ISSUES

A. System Performance

List samples qualified based on the degradation of system performance during sample analysis:

Sample ID	Comments	Actions
=====	=====	=====
_____	_____	_____
_____	_____	_____
_____	_____	_____
_____	_____	_____

Action:

Use professional judgment to qualify the data if it is determined that system performance has degraded during sample analyses. Inform the Contract Laboratory Program COR any action as a result of degradation of system performance which significantly affected the data.

B. Overall Assessment of Data

List samples qualified based on other issues:

Sample ID	Comments	Actions
=====	=====	=====
_____	_____	_____
No other issues that required the need to qualify the data. Results are valid and can be used for decision purposes.		
_____	_____	_____

Action:

1. Use professional judgment to determine if there is any need to qualify data which were not qualified based on the Quality Control (QC) criteria previously discussed.
2. Write a brief narrative to give the user an indication of the analytical limitations of the data. Inform the Contract Laboratory COR the action, any inconsistency of the data with the Sample Delivery Group (SDG) Narrative. If sufficient information on the intended use and required quality of the data is available, the reviewer should include their assessment of the usability of the data within the given context. This may be used as part of a formal Data Quality Assessment (DQA).
3. Sometimes, due to dilutions, re-analysis or SIM/Scan runs are being performed, there will be multiple results for a single analyte from a single sample. The following criteria and professional judgment are used to determine which result should be reported:
 - The analysis with the lower CRQL
 - The analysis with the better QC results
 - The analysis with the higher results

EXECUTIVE NARRATIVE

SDG No: **JC19473** Laboratory: **Accutest, Florida**
Analysis: **SW846-8015C** Number of Samples: **3**
Location: **BMSMC, Building 5 Area**
Humacao, PR

SUMMARY: Three (3) samples were analyzed for the low molecular weight alcohols (LMWAs) list following method SW846-8015C. The sample results were assessed according to USEPA data validation guidance documents in the following order of precedence: "Test Methods for Evaluating Solid Waste, Physical/Chemical Methods SW-846 (Final Update III, December 1996)," specifically for Methods 8000/8015C are utilized. The QC criteria and data validation actions listed on the data review worksheets are from the primary guidance document, unless otherwise noted.

Results are valid and can be used for decision making purposes.

Critical issues: **None**
Major: **None**
Minor: **None**

Critical findings: **None**
Major findings: **None**
Minor findings:

1. All samples properly preserved except for the followings: Sample JC19473-3 was received in a bulk container but was not preserved within 48 hours of sampling. Analysis performed at Accutest Laboratories, Orlando FL. Results qualified as (UJ) in affected samples.
2. Surrogate recovery outside laboratory control limit in sample JC19473-3. No action taken, % recovery within generally acceptable control limits. Target analytes not detected in the sample.
3. MSR recovery for several analytes outside laboratory control limits. No action taken, MS/MSD applies only the unspiked sample. Unspiked sample from another data package, was used for QC purposes only.
4. Blank spike recovery for n-Butyl Alcohol outside the laboratory control limit. No action taken, blank spike recovery within generally acceptable control limits.

COMMENTS: Results are valid and can be used for decision making purposes.

Reviewers Name: **Rafael Infante**
Chemist License 1888

Signature:



Date: **May 23, 2016**

SAMPLE ORGANIC DATA SAMPLE SUMMARY

Sample ID: JC19473-1

Sample location: BMSMC Building 5 Area

Sampling date: 4/30/2016

Matrix: AQ - Equipment Blank

METHOD: 8015C

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Ethanol	5.0	mg/l	1.0	-	U	Yes
Isobutyl Alcohol	5.0	mg/l	1.0	-	U	Yes
Isopropyl Alcohol	5.0	mg/l	1.0	-	U	Yes
n-Propyl Alcohol	5.0	mg/l	1.0	-	U	Yes
n-Butyl Alcohol	5.0	mg/l	1.0	-	U	Yes
Methanol	5.0	mg/l	1.0	-	U	Yes

Sample ID: JC19473-2A

Sample location: BMSMC Building 5 Area

Sampling date: 5/2/2016

Matrix: Groundwater

METHOD: 8015C

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Ethanol	5.0	mg/l	1.0	-	U	Yes
Isobutyl Alcohol	5.0	mg/l	1.0	-	U	Yes
Isopropyl Alcohol	5.0	mg/l	1.0	-	U	Yes
n-Propyl Alcohol	5.0	mg/l	1.0	-	U	Yes
n-Butyl Alcohol	5.0	mg/l	1.0	-	U	Yes
Methanol	5.0	mg/l	1.0	-	U	Yes

Sample ID: JC19473-3

Sample location: BMSMC Building 5 Area

Sampling date: 5/2/2016

Matrix: Soil

METHOD: 8015C

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Ethanol	14	mg/kg	1.0	-	UJ	Yes
Isobutyl Alcohol	14	mg/kg	1.0	-	UJ	Yes
Isopropyl Alcohol	14	mg/kg	1.0	-	UJ	Yes
n-Propyl Alcohol	14	mg/kg	1.0	-	UJ	Yes
n-Butyl Alcohol	14	mg/kg	1.0	-	UJ	Yes
Methanol	14.0	mg/kg	1.0	-	UJ	Yes

DATA REVIEW WORKSHEETS

Project Number: JC19473
 Date: 04/30-05/02/2016
 Shipping Date: 05/02/2016
 EPA Region: 2

REVIEW OF VOLATILE ORGANIC PACKAGE

The following guidelines for evaluating volatile organics were created to delineate required validation actions. This document will assist the reviewer in using professional judgment to make more informed decision and in better serving the needs of the data users. The sample results were assessed according to USEPA data validation guidance documents in the following order of precedence: "Test Methods for Evaluating Solid Waste, Physical/Chemical Methods SW-846 (Final Update III, December 1996)," specifically for Methods 8000/8015C are utilized. The QC criteria and data validation actions listed on the data review worksheets are from the primary guidance document, unless otherwise noted.

The hardcopied (laboratory name) Accutest data package received has been reviewed and the quality control and performance data summarized. The modified data review for VOCs included:

Lab. Project/SDG No.: JC19473 Sample matrix: Soil/Groundwater
 No. of Samples: 3

Trip blank No.: -
 Field blank No.: -
 Equipment blank No.: JC19473-1
 Field duplicate No.: -

<input checked="" type="checkbox"/> Data Completeness	<input checked="" type="checkbox"/> Laboratory Control Spikes
<input checked="" type="checkbox"/> Holding Times	<input checked="" type="checkbox"/> Field Duplicates
<input type="checkbox"/> N/A GC/MS Tuning	<input checked="" type="checkbox"/> Calibrations
<input type="checkbox"/> N/A Internal Standard Performance	<input checked="" type="checkbox"/> Compound Identifications
<input checked="" type="checkbox"/> Blanks	<input checked="" type="checkbox"/> Compound Quantitation
<input checked="" type="checkbox"/> Surrogate Recoveries	<input checked="" type="checkbox"/> Quantitation Limits
<input checked="" type="checkbox"/> Matrix Spike/Matrix Spike Duplicate	

Overall Comments: Low molecular weight alcohols by SW-846_8015C

Definition of Qualifiers:

J- Estimated results
 U- Compound not detected
 R- Rejected data
 UJ- Estimated nondetect

Reviewer: Rafael Dufur
 Date: May 23, 2016

DATA COMPLETENESS

DATE RECEIVED

Date		Page	
		1	
		2	
		3	
		4	
		5	
		6	
		7	
		8	
		9	
		10	
		11	
		12	
		13	
		14	
		15	
		16	
		17	
		18	
		19	
		20	
		21	
		22	
		23	
		24	
		25	
		26	
		27	
		28	
		29	
		30	
		31	
		32	
		33	
		34	
		35	
		36	
		37	
		38	
		39	
		40	
		41	
		42	
		43	
		44	
		45	
		46	
		47	
		48	
		49	
		50	
		51	
		52	
		53	
		54	
		55	
		56	
		57	
		58	
		59	
		60	
		61	
		62	
		63	
		64	
		65	
		66	
		67	
		68	
		69	
		70	
		71	
		72	
		73	
		74	
		75	
		76	
		77	
		78	
		79	
		80	
		81	
		82	
		83	
		84	
		85	
		86	
		87	
		88	
		89	
		90	
		91	
		92	
		93	
		94	
		95	
		96	
		97	
		98	
		99	
		100	

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below _____

HOLDING TIMES

The objective of this parameter is to ascertain the validity of the results based on the holding time of the sample from time of collection to the time of analysis.

Complete table for all samples and note the analysis and/or preservation not within criteria

SAMPLE ID	DATE SAMPLED	DATE ANALYZED	pH	ACTION
All samples analyzed within the recommended method holding time. All samples properly preserved except for the followings: Sample JC19473-3 was received in a bulk container but was not preserved within 48 hours of sampling. Analysis performed at Accutest Laboratories, Orlando FL. Results qualified as (UJ) in affected samples.				

Criteria

Aqueous samples – 14 days from sample collection for preserved samples ($\text{pH} \leq 2$, 4°C), no air bubbles.

Aqueous samples – 7 days from sample collection for unpreserved samples, 4°C , no air bubbles.

Soil samples- 7 days from sample collection.

Cooler temperature (Criteria: $4 \pm 2^{\circ}\text{C}$): 3.6°C

Actions

If the VOCs vial(s) have air bubbles, estimate positive results (J) and reject nondetects (R).

If the % solids of soil samples is 10-50%, estimate positive results (J) and nondetects (UJ)

If the % solid of soil samples is $< 10\%$, estimate positive results (J) and reject nondetects (R).

If holding times are exceeded but < 14 days beyond criteria, estimate positive results (J) and nondetects (UJ).

If holding times are exceeded but < 28 days beyond criteria, estimate positive results (J) and reject nondetects (R).

If holding times are grossly exceeded (> 28 days beyond criteria), reject all results (R).

If samples were not iced or if the ice were melted ($> 10^{\circ}\text{C}$), estimate positive results (J) and nondetects (UJ).

DATA REVIEW WORKSHEETS

All criteria were met __N/A__
Criteria were not met see below _____

GC/MS TUNING

The assessment of the tuning results is to determine if the sample instrumentation is within the standard tuning QC limits

__N/A__ The BFB performance results were reviewed and found to be within the specified criteria.

__N/A__ BFB tuning was performed for every 12 hours of sample analysis.

If no, use professional judgment to determine whether the associated data should be accepted, qualified or rejected.

List the samples affected: _____

If mass calibration is in error, all associated data are rejected.

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below _____

CALIBRATION VERIFICATION

Compliance requirements for satisfactory instrument calibration are established to ensure that the instrument is capable of producing and maintaining acceptable quantitative data.

Date of initial calibration: 05/02/16
 Dates of continuing calibration: 05/02/16 (initial); 05/05/16; 05/06/16
 Dates of final calibration verification: 05/02/16; 05/05/16; 05/06/16
 Instrument ID number: VOA5
 Matrix/Level: Aqueous/low

DATE	LAB FILE ID#	CRITERIA OUT RFs, %RSD, %D, r	COMPOUND	SAMPLES AFFECTED

Note: Initial, continuing, and final calibration verifications meets method specific criteria.

Criteria

All RFs must be > 0.05 regardless of method requirements for SPCC.

All %RSD must be $\leq 15\%$ regardless of method requirements for CCC.

All %Ds must be $\leq 20\%$ regardless of method requirements for CCC.

It should be noted that Region 2 SOP HW-24 does not specify criterion for the curve correlation coefficient (r). A limit for r of ≥ 0.995 has therefore been utilized as professional judgment.

Actions

If any compound has an initial RF or a continuing RF of < 0.05 , estimate positive results (J) and reject nondetects (R), regardless of method requirements.

If any compound has a %RSD $> 15\%$, estimate positive results (J) and use professional judgment to qualify nondetects.

If any compound has a %RSD $> 90\%$, estimate positive results (J) and reject nondetects (R).

If any compound has a % D $> 20\%$, estimate positive results (J) and reject nondetects (R).

If any compound has a % D $> 20\%$, estimate positive results (J) and nondetects (UJ).

If any compound has a % D $> 90\%$, estimate positive results (J) and reject nondetects (R).

If any compound has r < 0.995 , estimate positive results and nondetects.

A separate worksheet should be filled for each initial curve

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met _____
 and/or see below _____

V A. BLANK ANALYSIS RESULTS (Sections 1 & 2)

The assessment of the blank analysis results is to determine the existence and magnitude of contamination problems. The criteria for evaluation of blanks apply only to blanks associated with the samples, including trip, equipment, and laboratory blanks. If problems with any blanks exist, all data associated with the case must be carefully evaluated to determine whether or not there is an inherent variability in the data for the case, or if the problem is an isolated occurrence not affecting other data.

List the contamination in the blanks below. High and low levels blanks must be treated separately.

Laboratory blanks

DATE ANALYZED	LAB ID	LEVEL/MATRIX	COMPOUND	CONCENTRATION UNITS
_____	_____	_____	_____	_____
_____	_____	_____	_____	_____
All method blank meets method specific criteria				
_____	_____	_____	_____	_____
_____	_____	_____	_____	_____
_____	_____	_____	_____	_____
_____	_____	_____	_____	_____

Field/Equipment/Trip blank

DATE ANALYZED	LAB ID	LEVEL/MATRIX	COMPOUND	CONCENTRATION UNITS
No field/trip blanks included in this data package. No target analytes detected in the equipment blank.				
_____	_____	_____	_____	_____
_____	_____	_____	_____	_____
_____	_____	_____	_____	_____
_____	_____	_____	_____	_____
_____	_____	_____	_____	_____
_____	_____	_____	_____	_____
_____	_____	_____	_____	_____
_____	_____	_____	_____	_____
_____	_____	_____	_____	_____

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

V B. BLANK ANALYSIS RESULTS (Section 3)

Blank Actions

Action Levels (ALs) should be based upon the highest concentration of contaminant determined in any blank. Do not qualify any blank with another blank. The ALs for samples which have been diluted should be corrected for the sample dilution factor and/or % moisture, where applicable. No positive sample results should be reported unless the concentration of the compound in the samples exceeds the ALs:

ALs = 10x the amount of common contaminants (methylene chloride, acetone, 2-butanone, and toluene)

ALs = 5x for any other compounds

Specific actions are as follows:

If the concentration is < sample quantitation limit (SQL) and \leq AL, report the compound as not detected (U) at the SQL.

If the concentration is \geq SQL but \leq AL, report the compound as not detected (U) at the reported concentration.

If the concentration is \geq SQL and $>$ AL, report the concentration unqualified.

Notes:

High and low level blanks must be treated separately

Compounds qualified "U" for blank contamination are still considered "hits" when qualifying for calibration criteria.

CONTAMINATION SOURCE/LEVEL	COMPOUND	CONC/UNITS	AL/UNITS	SQL	AFFECTED SAMPLES

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

SURROGATE SPIKE RECOVERIES

Laboratory performance of individual samples is established by evaluation of surrogate spike recoveries. All samples are spiked with surrogate compounds prior to sample analysis. The accuracy of the analysis is measured by the surrogate percent recovery. Since the effects of the sample matrix are frequently outside the control of the laboratory and may present relatively unique problems, the validation of data is frequently subjective and demands analytical experience and professional judgment.

List the percent recoveries (%Rs) which do not meet the criteria for surrogate recovery.

Matrix: solid/aqueous

SAMPLE ID	SURROGATE COMPOUND				ACTION
	Hexanol	DBFM	TOL-d8	BFB	

 All surrogate recoveries within laboratory control limits except for the following:

<u> JC19473-3 </u>	<u> 130 % </u>	<u> No action </u>
----------------------	-------------------	----------------------

Note: No action taken, % recovery within generally acceptable control limits. Target analytes not detected in the sample

QC Limits* (Aqueous)

 LL to UL 73 to 123 to to to

QC Limits* (Solid-Low)

 LL to UL 69 to 121 to to to

QC Limits* (Solid-Med)

 LL to UL to to to to

1,2-DCA = 1,2-Dichloromethane-d4

TOL-d8 = Toluene-d8

DBFM = Dibromofluoromethane

BFB = Bromofluorobenzene

* QC limits are laboratory in-house performance criteria, LL = lower limit, UL = upper limit.

* If QC limits are not available, use limits of 80 – 120 % for aqueous and 70 – 130 % for solid samples.

Actions:

QUALITY	%R < 10%	%R = 10% - LL	%R > UL
Positive results	J	J	J
Nondetects results	R	UJ	Accept

Surrogate action should be applied:

If one or more surrogate in the VOC fraction is out of specification, but has a recovery of > 10%.

If any one surrogate in a fraction shows < 10 % recovery.

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

VII. A MATRIX SPIKE/MATRIX SPIKE DUPLICATE (MS/MSD)

This data is generated to determine long term precision and accuracy in the analytical method for various matrices. This data alone cannot be used to evaluate the precision and accuracy of individual samples. If any % R in the MS or MSD falls outside the designated range, the reviewer should determine if there are matrix effects, i.e. LCS data are within the QC limits but MS/MSD data are outside QC limit.

1. MS/MSD Recoveries and Precision Criteria

The laboratory should use one MS and a duplicate analysis of an unspiked field sample if target analytes are expected in the sample. If target analytes are not expected, MS/MSD should be analyzed.

List the %Rs, RPD of the compounds which do not meet the criteria.

Sample ID: JC19423-2AMS/-MSD Matrix/Level: Aqueous
 Sample ID: JC19423-1MS/-1MSD Matrix/Level: Soil

MS OR MSD	COMPOUND	% R	RPD	QC LIMITS	ACTION
-----------	----------	-----	-----	-----------	--------

MS/MSD%_recoveries_and_RPD_within_laboratory_control_limits_except_for_the_followings:_

The QC reported here applies to the following samples: Method: SW846 8015C MOD
 JC19473-3

Compound	JC19423-1 mg/kg	Q	Spike mg/kg	MS mg/kg	MS %	Spike mg/kg	MSD mg/kg	MSD %	RPD	Limits Rec/RPD
Ethanol	ND		282	298	106	282	331	118*	10	80-117/13
Isopropyl Alcohol	ND		282	303	108	282	340	121*	12	75-116/15
n-Propyl Alcohol	ND		282	301	107	282	337	120*	11	78-116/13
n-Butyl Alcohol	ND		282	303	108	282	339	120*	11	74-115/13

Surrogate Recoveries	MS	MSD	JC19423-1	Limits
Hexanol	102%	114%	112%	69-121%

(a) Sample was received in a bulk container but was not preserved within 48 hours of sampling.

- * QC limits are laboratory in-house performance criteria, LL = lower limit, UL = upper limit.
- * If QC limits are not available, use limits of 70 – 130 %.

Actions:

QUALITY	%R < LL	%R > UL
Positive results	J	J
Nondetects results	R	Accept

DATA REVIEW WORKSHEETS

All criteria were met ☒ X
 Criteria were not met
 and/or see below _____

MS/MSD criteria apply only to the unspiked sample, its dilutions, and the associated MS/MSD samples:

If the % R for the affected compounds were < LL (or 70 %), qualify positive results (J) and nondetects (UJ).

If the % R for the affected compounds were > UL (or 130 %), only qualify positive results (J).

If 25 % or more of all MS/MSD %R were < LL (or 70 %) or if two or more MS/MSD %Rs were < 10%, qualify all positive results (J) and reject nondetects (R).

Note: No action taken, MS/MSD applies only the unspiked sample. Unspiked sample from another data package, used for QC purposes only

VII. B MATRIX SPIKE/MATRIX SPIKE DUPLICATE

MS/MSD – Unspiked Compounds

It should be noted that Region 2 SOP HW-24 does not specify a MS/MSD criteria for the unspiked compounds in the sample. A %RSD of < 50% has therefore been utilized as professional judgment.

If all target analytes were spiked in the MS/MSD, this review element is not applicable.

List the %RSD of the compounds which do not meet the criteria.

Sample ID: _____ - _____ Matrix/Level/Unit: _____ - _____

COMPOUND	SAMPLE CONC.	MS CONC.	MSD CONC.	% RSD	ACTION

Actions:

* If the % RSD > 50, qualify the positive result in the unspiked samples as estimated (J).

* If the % RSD is not calculated (NC) due to nondetected value, use professional judgment to qualify the data.

A separate worksheet should be used for each MS/MSD pair.

DATA REVIEW WORKSHEETS

All criteria were met X
Criteria were not met
and/or see below

VIII. LABORATORY CONTROL SAMPLE (LCS) ANALYSIS

This data is generated to determine accuracy of the analytical method for various matrices.

1. LCS Recoveries Criteria

Where LCS spiked with the same analyte at the same concentrations as the MS/MSD? **Yes**
or No. If no make note in data review memo.

List the %R of compounds which do not meet the criteria

LCS ID	COMPOUND	% R	QC LIMIT
Recoveries within laboratory control limits except for the following:			
GXY2775-BS_(Solid)	n-Butyl Alcohol	116 %	74 - 115

Note: No action taken, blank spike recovery within generally acceptable control limits.

- * QC limits are laboratory in-house performance criteria, LL = lower limit, UL = upper limit.
- * If QC limits are not available, use limits of 70 – 130 %.

Actions:

QUALITY	%R < LL	%R > UL
Positive results	J	J
Nondetects results	R	Accept

All analytes in the associated sample results are qualified for the following criteria.

If 25 % of the LCS recoveries were < LL (or 70 %), qualify all positive results (j) and reject nondetects (R).

If two or more LCS were below 10 %, qualify all positive results as (J) and reject nondetects (R).

2. Frequency Criteria:

Where LCS analyzed at the required frequency and for each matrix? Yes or No.

If no, the data may be affected. Use professional judgment to determine the severity of the effect and qualify data accordingly. Discuss any actions below and list the samples affected.

DATA REVIEW WORKSHEETS

All criteria were met N/A
 Criteria were not met
 and/or see below

IX. FIELD/LABORATORY DUPLICATE PRECISION

Sample IDs: - Matrix: -

Field/laboratory duplicates samples may be taken and analyzed as an indication of overall precision. These analyses measure both field and lab precision; therefore, the results may have more variability than laboratory duplicates which only laboratory performance. It is also expected that soil duplicate results will have a greater variance than water matrices due to difficulties associated with collecting identical field duplicate samples.

The project QAPP should be reviewed for project-specific information.

Suggested criteria: RPD \pm 30% for aqueous samples, RPD \pm 50 % for solid samples. If both samples and duplicate are <5 SQL, the RPD criteria is doubled.

COMPOUND	SQL	SAMPLE CONC.	DUPLICATE CONC.	RPD	ACTION
No field/laboratory duplicate analyzed with this data package. MS/MSD % recoveries RPD used to assess precision. RPD within laboratory and generally acceptable control limits.					

Actions:

Qualify as estimated positive results (J) and nondetects (UJ) for the compound that exceeded the above criteria. For organics, only the sample and duplicate will be qualified.

If an RPD cannot be calculated because one or both of the sample results is not detected, the following actions apply:

If one sample result is not detected and the other is greater than 5x the SQL qualify (J/UJ).

If one sample value is not detected and the other is greater than 5x the SQL and the SQLs for the sample and duplicate are significantly different, use professional judgment to determine if qualification is appropriate.

If one sample value is not detected and the other is less than 5x, use professional judgment to determine if qualification is appropriate.

If both sample and duplicate results are not detected, no action is needed.

The assessment of the internal standard (IS) parameter is used to assist the data reviewer in determining the condition of the analytical instrumentation.

* Area of +100% or -50% of the IS area in the associated calibration standard.
* Retention time (RT) within 30 seconds of the IS area in the associated calibration standard.

[illegible]

1. IS actions should be applied to the compound quantitated with the out-of-control ISs

QUALITY	IS AREA < -25%	IS AREA = -25 % TO - 50%	IS AREA > + 100%
Positive results	J	J	J
Nondetected results	R	UJ	ACCEPT

- 13

DATA REVIEW WORKSHEETS

All criteria were met X
Criteria were not met
and/or see below _____

XII. SAMPLE QUANTITATION

The sample quantitation evaluation is to verify laboratory quantitation results. In the space below, please show a minimum of one sample calculation:

JC19473-1

Hexanol

RF = 5304

$$[] = (482907)/(5304)$$

$$= 91.04 \text{ ppm OK}$$

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

XII. QUANTITATION LIMITS

A. Dilution performed

SAMPLE ID	DILUTION FACTOR	REASON FOR DILUTION

B. Percent Solids

List samples which have ≤ 50 % solids

Actions:

If the % solids of a soil sample is 10-50%, estimate positive results (J) and nondetects (UJ)

If the % solids of a soil sample is < 10%, estimate positive results (J) and reject nondetects (R)

EXECUTIVE NARRATIVE

SDG No: **JC19473** Laboratory: **Accutest, New Jersey**
Analysis: **SW846-8081B** Number of Samples: **3**

Location: **BMSMC, Building 5 Area**
Humacao, PR

SUMMARY: Three (3) samples were analyzed for selected pesticides following method SW846-8081B. The sample results were assessed according to USEPA data validation guidance documents in the following order of precedence *Hazardous Waste Support Section SOP No. HW-36A, Revision 0, June, 2015. SOM02.2. Pesticide Data Validation*. The QC criteria and data validation actions listed on the data review worksheets are from the primary guidance document, unless otherwise noted.

Results are valid and can be used for decision making purposes.

Critical issues: **None**
Major issues: **None**
Minor issues: **None**

Critical findings: **None**
Major findings: **None**
Minor findings: 1. **MS/MSD % recoveries and RPD outside laboratory control limits for several analytes. No action taken, MS/MSD results apply only to the unspiked sample.**
2. **No information for florisol cartridge performance check included in data package. Florisol cartridge used for sample extraction/clean-up. No qualification of the data performed, professional judgment.**

COMMENTS: Results are valid and can be used for decision making purposes.

Reviewers Name: **Rafael Infante**
Chemist License 1888

Signature:



Date: **May 23, 2016**

SAMPLE ORGANIC DATA SAMPLE SUMMARY

Sample ID: JC19473-1

Sample location: BMSMC Building 5 Area

Sampling date: 30-Apr-16

Matrix: AQ - Equipment Blank

METHOD: 8081B

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Aldrin	0.0067	ug/L	1	-	U	Yes
alpha-BHC	0.0067	ug/L	1	-	U	Yes
beta-BHC	0.0067	ug/L	1	-	U	Yes
delta-BHC	0.0067	ug/L	1	-	U	Yes
gamma-BHC (Lindane)	0.0067	ug/L	1	-	U	Yes
alpha-Chlordane	0.0067	ug/L	1	-	U	Yes
gamma-Chlordane	0.0067	ug/L	1	-	U	Yes
Dieldrin	0.0067	ug/L	1	-	U	Yes
4,4'-DDD	0.0067	ug/L	1	-	U	Yes
4,4'-DDE	0.0067	ug/L	1	-	U	Yes
4,4'-DDT	0.0067	ug/L	1	-	UJ	Yes
Endrin	0.0067	ug/L	1	-	U	Yes
Endosulfan sulfate	0.0067	ug/L	1	-	U	Yes
Endrin aldehyde	0.0067	ug/L	1	-	U	Yes
Endrin ketone	0.0067	ug/L	1	-	U	Yes
Endosulfan-I	0.0067	ug/L	1	-	U	Yes
Endosulfan-II	0.0067	ug/L	1	-	U	Yes
Heptachlor	0.0067	ug/L	1	-	U	Yes
Heptachlor epoxide	0.0067	ug/L	1	-	U	Yes
Methoxychlor	0.013	ug/L	1	-	UJ	Yes
Toxaphene	0.17	ug/L	1	-	U	Yes

Sample ID: JC19473-2
Sample location: BMSMC Building 5 Area
Sampling date: 2-May-16
Matrix: Groundwater

METHOD: 8081B

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Aldrin	0.0071	ug/l	1	-	U	Yes
alpha-BHC	0.0071	ug/l	1	-	U	Yes
beta-BHC	0.0071	ug/l	1	-	U	Yes
delta-BHC	0.0071	ug/l	1	-	U	Yes
gamma-BHC (Lindane)	0.0071	ug/l	1	-	U	Yes
alpha-Chlordane	0.0071	ug/l	1	-	U	Yes
gamma-Chlordane	0.0071	ug/l	1	-	U	Yes
Dieldrin	0.0071	ug/l	1	-	U	Yes
4,4'-DDD	0.0071	ug/l	1	-	U	Yes
4,4'-DDE	0.0071	ug/l	1	-	U	Yes
4,4'-DDT	0.0071	ug/l	1	-	U	Yes
Endrin	0.0071	ug/l	1	-	U	Yes
Endosulfan sulfate	0.0071	ug/l	1	-	U	Yes
Endrin aldehyde	0.0071	ug/l	1	-	U	Yes
Endrin ketone	0.0071	ug/l	1	-	U	Yes
Endosulfan-I	0.0071	ug/l	1	-	U	Yes
Endosulfan-II	0.0071	ug/l	1	-	U	Yes
Heptachlor	0.0071	ug/l	1	-	U	Yes
Heptachlor epoxide	0.0071	ug/l	1	-	U	Yes
Methoxychlor	0.014	ug/l	1	-	U	Yes
Toxaphene	0.18	ug/l	1	-	U	Yes

Sample ID: JC19473-3
Sample location: BMSMC Building 5 Area
Sampling date: 2-May-16
Matrix: Soil

METHOD: 8081B

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Aldrin	0.86	ug/kg	1	-	U	Yes
alpha-BHC	0.86	ug/kg	1	-	U	Yes
beta-BHC	0.86	ug/kg	1	-	U	Yes
delta-BHC	0.86	ug/kg	1	-	U	Yes
gamma-BHC (Lindane)	0.86	ug/kg	1	-	U	Yes
alpha-Chlordane	0.86	ug/kg	1	-	U	Yes
gamma-Chlordane	0.86	ug/kg	1	-	U	Yes
Dieldrin	0.86	ug/kg	1	-	U	Yes
4,4'-DDD	0.86	ug/kg	1	-	U	Yes
4,4'-DDE	0.86	ug/kg	1	-	U	Yes
4,4'-DDT	0.86	ug/kg	1	-	U	Yes
Endrin	0.86	ug/kg	1	-	U	Yes
Endosulfan sulfate	0.86	ug/kg	1	-	U	Yes
Endrin aldehyde	0.86	ug/kg	1	-	U	Yes
Endosulfan-I	0.86	ug/kg	1	-	U	Yes
Endosulfan-II	0.86	ug/kg	1	-	U	Yes
Heptachlor	0.86	ug/kg	1	-	U	Yes
Heptachlor epoxide	0.86	ug/kg	1	-	U	Yes
Methoxychlor	1.7	ug/kg	1	-	U	Yes
Endrin ketone	0.86	ug/kg	1	-	U	Yes
Toxaphene	21	ug/kg	1	-	U	Yes

DATA REVIEW WORKSHEETS

Project/Case Number: JC19473
 Sampling Date: 04/30-05/02, 2016
 Shipping Date: May 02, 2016
 EPA Region No.: 2

REVIEW OF PESTICIDE ORGANIC PACKAGE

The following guidelines for evaluating volatile organics were created to delineate required validation actions. This document will assist the reviewer in using professional judgment to make more informed decision and in better serving the needs of the data users. The sample results were assessed according to USEPA data validation guidance documents in the following order of precedence *Hazardous Waste Support Section SOP No. HW-36A, Revision 0, June, 2015. SOM02.2. Pesticide Data Validation*. The QC criteria and data validation actions listed on the data review worksheets are from the primary guidance document, unless otherwise noted.

The hardcopied (laboratory name) Accutest data package received has been reviewed and the quality control and performance data summarized. The data review for VOCs included:

Lab. Project/SDG No.: JC19473 Sample matrix: SoilGroundwater
 No. of Samples: 3

Trip blank No.: -
 Field blank No.: -
 Equipment blank No.: JC19473-1
 Field duplicate No.: -
 Field spikes No.: -
 QC audit samples: -

<input checked="" type="checkbox"/> Data Completeness	<input checked="" type="checkbox"/> Laboratory Control Spikes
<input checked="" type="checkbox"/> Holding Times	<input checked="" type="checkbox"/> Field Duplicates
<input type="checkbox"/> GC/MS Tuning	<input checked="" type="checkbox"/> Calibrations
<input checked="" type="checkbox"/> Internal Standard Performance	<input checked="" type="checkbox"/> Compound Identifications
<input checked="" type="checkbox"/> Blanks	<input checked="" type="checkbox"/> Compound Quantitation
<input checked="" type="checkbox"/> Surrogate Recoveries	<input checked="" type="checkbox"/> Quantitation Limits
<input checked="" type="checkbox"/> Matrix Spike/Matrix Spike Duplicate	

Overall Comments: TCL_pesticides_list_by_SW846-8081B

Definition of Qualifiers:

J- Estimated results
 U- Compound not detected
 R- Rejected data
 UJ- Estimated nondetect

Reviewer: Rafael Delacruz
 Date: May 23, 2016

DATA COMPLETENESS

DATE RECEIVED

This image shows a single sheet of white paper with horizontal blue or grey ruling lines. A diagonal line runs from the top left corner towards the bottom right, creating a narrow margin on the left side. The paper appears to be a template for writing or drawing.

DATA REVIEW WORKSHEETS

All criteria were met X
Criteria were not met
and/or see below

HOLDING TIMES

The objective of this parameter is to ascertain the validity of the results based on the holding time of the sample from time of collection to the time of analysis.

Complete table for all samples and note the analysis and/or preservation not within criteria

SAMPLE ID	DATE SAMPLED	DATE EXTRACTED/ANALYZED	ACTION
Samples properly preserved.			

Preservatives: All samples extracted and analyzed within the required criteria.

Criteria

Aqueous samples - seven (7) days from sample collection for extraction; 40 days from sample collection for analysis.

Non-aqueous samples – fourteen (14) days from sample collection for extraction; 40 days from sample collection for analysis.

Cooler temperature (Criteria: 4 ± 2 °C): 3.6°C - OK

Actions

Qualify aqueous sample results using preservation and technical holding time information as follows:

- If there is no evidence that the samples were properly preserved ($T = 4^{\circ}\text{C} \pm 2^{\circ}\text{C}$), and the samples were extracted or analyzed within the technical holding times, qualify detects as estimated (J) and non-detects as estimated (UJ).
- If there is no evidence that the samples were properly preserved ($T = 4^{\circ}\text{C} \pm 2^{\circ}\text{C}$), and the samples were extracted or analyzed outside the technical holding times, qualify detects as estimated (J) and non-detects as estimated (UJ).
- If the samples were properly preserved, and were extracted and analyzed within the technical holding times, no qualification of the data is necessary.
- If the samples were properly preserved, and were extracted or analyzed outside the technical holding times, qualify detects as estimated (J) and non-detects as estimated (UJ). Note in the Data Review Narrative that holding times were exceeded and the effect of exceeding the holding time on the resulting data.

DATA REVIEW WORKSHEETS

- e. Use professional judgment to qualify samples whose temperature upon receipt at the laboratory is either below 2 degrees centigrade or above 6 degrees centigrade.
- f. If technical holding times are grossly exceeded, use professional judgment to qualify the data.

Qualify non-aqueous sample results using preservation and technical holding time information as follows:

- a. If there is no evidence that the samples were properly preserved ($T = 4^{\circ}\text{C} \pm 2^{\circ}\text{C}$), and the samples were extracted or analyzed within the technical holding time, qualify detects as estimated (J) and non-detects as estimated (UJ).
- b. If there is no evidence that the samples were properly preserved ($T = 4^{\circ}\text{C} \pm 2^{\circ}\text{C}$), and the samples were extracted or analyzed outside the technical holding time, qualify detects as estimated (J) and non-detects as estimated (UJ).
- c. If the samples were properly preserved, and were extracted and analyzed within the technical holding time, no qualification of the data is necessary.
- d. If the samples were properly preserved, and were extracted or analyzed outside the technical holding time, qualify detects as estimated (J) and non-detects as estimated (UJ). Note in the Data Review Narrative that holding times were exceeded and the effect of exceeding the holding time on the resulting data.
- e. Use professional judgment to qualify samples whose temperature upon receipt at the laboratory is either below 2 degrees centigrade or above 6 degrees centigrade.
- f. If technical holding times are grossly exceeded, use professional judgment to qualify the data.

All criteria were met X
Criteria were not met see below

**GAS CHROMATOGRAPH WITH ELECTRON CAPTURE DETECTOR (GC/ECD) INSTRUMENT
PERFORMANCE CHECK (SECTIONS 1 TO 5)**

1. Resolution Check Mixture

Criteria

Is the resolution between two adjacent peaks in the Resolution Check Mixture C greater than or equal to 80.0% for all analytes for the primary column and greater than or equal to 50.0% for the confirmation column? Yes? or No?

Is the resolution between two adjacent peaks in the Resolution Check Mixture (A and B) greater than or equal to 60.0%? Yes? or No?

Note: If resolution criteria are not met, the quantitative results may not be accurate due to inadequate resolution. Qualitative identifications may also be questionable if coelution exists.

Action

- a. Qualify detects for target compounds that were not adequately resolved as tentatively identified (NJ).
- b. Qualify non-detected compounds as unusable (R).

2. Performance Evaluation Mixture (PEM) Resolution Criteria

Criteria

Is PEM analysis performed at the required frequency (at the end of each pesticide initial calibration sequence and every 12 hours)? Yes? or No?

Action

- a. If PEM is not performed at the required frequency, qualify all associated sample and blank results as unusable (R).

Criteria

Is PEM % Resolution < 90%? Yes? or No?

Action

- a. a. Qualify detects for target compounds that were not adequately resolved as tentatively identified (NJ).
- b. Qualify non-detected compounds as unusable (R).

DATA REVIEW WORKSHEETS

All criteria were met X
Criteria were not met see below

3. PEM 4,4'-DDT Breakdown

Criteria

Is the PEM 4,4'-DDT % Breakdown >20.0% and 4,4'-DDT is detected? Yes? or No?

Action

a. Qualify detects for 4,4'-DDT; detects for 4,4'-DDD; and detects for 4,4'-DDE as estimated (J)

Criteria

Is the PEM 4,4'-DDT % Breakdown >20.0% and 4,4'-DDT is not detected Yes? or No?

Action

- a. Qualify non-detects for 4,4'- DDT as unusable (R)
- b. Qualify detects for 4,4'-DDD as tentatively identified (NJ)
- c. Qualify detects for 4,4'-DDE as tentatively identified (NJ)

4. PEM Endrin Breakdown

Criteria

Is the PEM Endrin % Breakdown >20.0% and Endrin is detected? Yes? or No?

Action

a. Qualify detects for Endrin; detects for Endrin aldehyde; and detects for Endrin ketone as estimated (J)

Criteria

Is the PEM Endrin % Breakdown >20.0% and Endrin is not detected Yes? or No?

Action

- a. Qualify non-detects for Endrin as unusable (R)
- b. Qualify detects for Endrin aldehyde as tentatively identified (NJ)
- c. Qualify detects for Endrin ketone as tentatively identified (NJ)

DATA REVIEW WORKSHEETS

All criteria were met X
Criteria were not met see below

5. Mid-point Individual Standard Mixture Resolution -

Criteria

Is the resolution between two adjacent peaks in the Resolution Check Mixture C greater than or equal to 80.0% for all analytes for the primary column and greater than or equal to 50.0% for the confirmation column?
Yes? or No?

Is the resolution between two adjacent peaks in the Resolution Check Mixture (A and B) greater than or equal to 90.0%?
Yes? or No?

Note: If resolution criteria are not met, the quantitative results may not be accurate due to inadequate resolution. Qualitative identifications may also be questionable if coelution exists.

Action

- a. Qualify detects for target compounds that were not adequately resolved as tentatively identified (NJ).
- b. Qualify non-detected compounds as unusable (R).

Criteria

Is mid-point individual standard mixture analysis performed at the required frequency (every 12 hours)?
Yes? or No?

Action

- a. If the mid-point individual standard mixture analysis is not performed at the required frequency, qualify all associated sample and blank results as unusable (R).

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below _____

CALIBRATION VERIFICATION

Compliance requirements for satisfactory instrument calibration are established to ensure that the instrument is capable of producing and maintaining acceptable quantitative data.

Date of initial calibration: 05/02/16 04/14/16
 Dates of initial calibration verification: 05/02/16 04/14/16
 Dates of continuing calibration: 05/06/16; 05/07/16 05/03/16; 05/04/16
 Dates of final calibration - -
 Instrument ID numbers: GC4G HP_1530A
 Matrix/Level: Aqueous/low Aqueous/low

DATE	LAB ID#	FILE	CRITERIA OUT RFs, %RSD, %D, r	COMPOUND	SAMPLES AFFECTED
Initial and initial calibration verification within the guidance document performance criteria. Continuing calibration % differences meet the performance criteria in at least one of the column. No final calibration verification analyzed. No action taken, professional judgment					

Criteria

Are a five point calibration curve delivered with concentration levels as shown in Table 3 of SOP HW-36A, Revision 0, June, 2015? Yes? or No?

Actions

If the standard concentrations listed in Table 3 are not used, use professional judgment to evaluate the effect on the data

Criteria

Are RT Windows calculated correctly? Yes? or No?

Action

Recalculate the windows and use the corrected values for all evaluations.

Criteria

Are the Percent Relative Standard Deviation (%RSD) of the CFs for each of the single component target compounds less than or equal to 20.0%, except for alpha-BHC and delta-BHC?

Yes? or No?

DATA REVIEW WORKSHEETS

Are the %RSD of the CFs for alpha-BHC and delta-BHC less than or equal to 25.0%. Yes? or No?

Is the %RSD of the CFs for each of the Toxaphene peaks must be < 30% when 5-point ICAL is performed? Yes? or No?

Is the %RSD of the CFs for the two surrogates (tetrachloro-m-xylene and decachlorobiphenyl) less than or equal to 30.0%. Yes? or No?

Action

- If the %RSD criteria are not met, qualify detects as estimated (J) and use professional judgment to qualify non-detected target compounds.
- If the %RSD criteria are within allowable limits, no qualification of the data is necessary

Continuing Calibration Checks

Criteria

Is the continuing calibration standard analyzed at the acceptable time intervals? Yes? or No?

Action

- If more than 14 hours has elapsed from the injection of the instrument blank that begins an analytical sequence (opening CCV) and the injection of either a PEM or mid-point concentration of the Individual Standard Mixtures (A and B) or (C), qualify all data as unusable (R).
- If more than 12 hours has elapsed from the injection of the instrument blank that begins an analytical sequence (opening CCV) and the injection of the last sample or blank that is part of the same analytical sequence, qualify all data as unusable (R).
- If more than 72 hours has elapsed from the injection of the sample with a Toxaphene detection and the Toxaphene Calibration Verification Standard (CS3), qualify all data as unusable (R).

Criteria

Is the Percent Difference (%D) within $\pm 25.0\%$ for the PEM sample? Yes? or No?

Action

- Qualify associated detects as estimated (J) and non-detects as estimated (UJ).

Criteria

For the Calibration Verification Standard (CS3); is the Percent Difference (%D) within $\pm 25.0\%$? Yes? or No?

Action

Qualify associated detects as estimated (J) and non-detects as estimated (UJ).

DATA REVIEW WORKSHEETS

All criteria were met X
Criteria were not met
and/or see below _____

Criteria

Is the PEM 4,4'-DDT % Breakdown >20.0% and 4,4'-DDT is detected? Yes? or No?

Action

- a. Qualify detects for 4,4'-DDT; detects for 4,4'-DDD; and detects for 4,4'-DDE as estimated (J)
- b. Non-detected associated compounds are not qualified

Criteria

Is the PEM 4,4'-DDT % Breakdown >20.0% and 4,4'-DDT is not detected Yes? or No?

Action

- a. Qualify non-detects for 4,4'- DDT as unusable (R)
- b. Qualify detects for 4,4'-DDD as tentatively identified (NJ)
- c. Qualify detects for 4,4'-DDE as tentatively identified (NJ)

Criteria

Is the PEM Endrin % Breakdown >20.0% and Endrin is detected? Yes? or No?

Action

- a. Qualify detects for Endrin; detects for Endrin aldehyde; and detects for Endrin ketone as estimated (J)
- b. Non-detected associated compounds are not qualified

Criteria

Is the PEM Endrin % Breakdown >20.0% and Endrin is not detected Yes? or No?

Action

- a. Qualify non-detects for Endrin as unusable (R)
- b. Qualify detects for Endrin aldehyde as tentatively identified (NJ)
- c. Qualify detects for Endrin ketone as tentatively identified (NJ)

A separate worksheet should be filled for each initial curve

DATA REVIEW WORKSHEETS

All criteria were met X
Criteria were not met
and/or see below _____

BLANK ANALYSIS RESULTS (Sections 1 & 2)

The assessment of the blank analysis results is to determine the existence and magnitude of contamination problems. The criteria for evaluation of blanks apply only to blanks associated with the samples, including trip, equipment, and laboratory blanks. If problems with any blanks exist, all data associated with the case must be carefully evaluated to determine whether or not there is an inherent variability in the data for the case, or if the problem is an isolated occurrence not affecting other data.

List the contamination in the blanks below. High and low levels blanks must be treated separately.

CRQL concentration N/A

Laboratory blanks

DATE ANALYZED	LAB ID	LEVEL/ MATRIX	COMPOUND	CONCENTRATION UNITS
---------------	--------	---------------	----------	---------------------

No target analytes detected in method blanks at a reporting limit of 0.01 and 0.001 ug/L.

Field/Equipment/Trip blank

DATE ANALYZED	LAB ID	LEVEL/ MATRIX	COMPOUND	CONCENTRATION UNITS
---------------	--------	---------------	----------	---------------------

No target analytes detected in the equipment blank. No field/trip blanks analyzed with this data package.

All criteria were met X
 Criteria were not met
 and/or see below

BLANK ANALYSIS RESULTS (Section 3)

Blank Actions

Action Levels (ALs) should be based upon the highest concentration of contaminant determined in any blank. Do not qualify any blank with another blank. The ALs for samples which have been diluted should be corrected for the sample dilution factor and/or % moisture, where applicable. No positive sample results should be reported unless the concentration of the compound in the samples exceeds the ALs:

The concentration of non-target compounds in all blanks must be less than or equal to 10 µg/L. The concentration of each target compound found in the method or field blanks must be less than its CRQL listed in the method.

Data concerning the field blanks are not evaluated as part of the CCS process. If field blanks are present, the data reviewer should evaluate this data in a similar fashion as the method blanks.

Specific actions are as follows:

Blank Actions for Pesticide Analyses

Blank Type	Blank Result	Sample Result	Action for Samples
Method, Sulfur Cleanup, Instrument, Field, TCLP/SPLP	Detects	Not detected	No qualification required
	< CRQL	< CRQL	Report CRQL value with a U
		≥ CRQL	No qualification required
	> CRQL	< CRQL	Report CRQL value with a U
		≥ CRQL and ≤ blank concentration	Report blank value for sample concentration with a U
		≥ CRQL and > blank concentration	No qualification required
	= CRQL	≤ CRQL	Report CRQL value with a U
		> CRQL	No qualification required
	Gross contamination	Detects	Report blank value for sample concentration with a U

• •

DATA REVIEW WORKSHEETS

All criteria were met ☒ X
 Criteria were not met
 and/or see below _____

SURROGATE SPIKE RECOVERIES

Laboratory performance of individual samples is established by evaluation of surrogate spike recoveries. All samples are spiked with surrogate compounds prior to sample analysis. The accuracy of the analysis is measured by the surrogate percent recovery. Since the effects of the sample matrix are frequently outside the control of the laboratory and may present relatively unique problems, the validation of data is frequently subjective and demands analytical experience and professional judgment.

List the percent recoveries (%Rs) which do not meet the criteria for surrogate recovery.

Matrix: Aqueous

Lab Sample ID	Lab File ID	S1 a	S1 b	S2 a	S2 b
JC19473-1	6G34643.D	57	58	54	50
JC19473-2	6G34644.D	77	78	103	84
OP93571-BS1	6G34622.D	81	80	85	93
OP93571-MB1	6G34621.D	84	83	76	83
OP93571-MS	6G34665.D	76	79	103	107
OP93571-MSD	6G34666.D	98	100	128* c	127* c

Surrogate Compounds	Recovery Limits
---------------------	-----------------

S1 = Tetrachloro-m-xylene	26-132%
---------------------------	---------

S2 = Decachlorobiphenyl	10-118%
-------------------------	---------

(a) Recovery from GC signal #1

(b) Recovery from GC signal #2

(c) Outside the QC limits.

Note: Surrogate recoveries within laboratory control limits in project samples. Outside control limit in QC sample (MSD). No action taken

DATA REVIEW WORKSHEETS

Matrix: Soil

Lab Sample ID	Lab File ID	S1 a	S1 b	S2 a	S2 b
JC19473-3	4G68007.D	110	107	139	108
OP93683-BS1	4G67964.D	109	88	131	97
OP93683-MB1	4G67963.D	104	87	127	94
OP93683-MS	4G67998.D	105	92	121	103
OP93683-MSD	4G67999.D	91	79	105	89

Surrogate Compounds Recovery Limits

S1 = Tetrachloro-m-xylene 24-136%

S2 = Decachlorobiphenyl 10-153%

(a) Recovery from GC signal #1

(b) Recovery from GC signal #2

Actions:

- For any surrogate recovery greater than 150%, qualify detected target compounds as biased high (J+).
- Do not qualify non-detected target compounds for surrogate recovery > 150 %.
- If both surrogate recoveries are greater than or equal to 30% and less than or equal to 150%, no qualification of the data is necessary.
- For any surrogate recovery greater than or equal to 10% and less than 30%, qualify detected target compounds as biased low (J-).
- For any surrogate recovery greater than or equal to 10% and less than 30%, qualify non-detected target compounds as approximated (UJ).
- If low surrogate recoveries are from sample dilution, professional judgment should be used to determine if the resulting data should be qualified. If sample dilution is not a factor:
 - Qualify detected target compounds as biased low (J-).
 - Qualify non-detected target compounds as unusable (R).
- If surrogate RTs in PEMs, Individual Standard Mixtures, samples, and blanks are outside of the RT Windows, the reviewer must use professional judgment to qualify data.
- If surrogate RTs are within RT windows, no qualification of the data is necessary.
- If the two surrogates were not added to all samples, MS/MSDs, standards, LCSs, and blanks, use professional judgment in qualifying data as missing surrogate analyte may not directly apply to target analytes.

DATA REVIEW WORKSHEETS

Summary Surrogate Actions for Pesticide Analyses

Criteria	Action*	
	Detected Target Compounds	Non-detected Target Compounds
%R > 150%	J+	No qualification
30% < %R < 150%	No qualification	
10% < %R < 30%	J-	UJ
%R < 10% (sample dilution not a factor)	J-	R
%R < 10% (sample dilution is a factor)	Use professional judgment	
RT out of RT window	Use professional judgment	
RT within RT window	No qualification	

- * Use professional judgment in qualifying data, as surrogate recovery problems may not directly apply to target analytes.

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

MATRIX SPIKE/MATRIX SPIKE DUPLICATE (MS/MSD)

This data is generated to determine long term precision and accuracy in the analytical method for various matrices. This data alone cannot be used to evaluate the precision and accuracy of individual samples. If any % R in the MS or MSD falls outside the designated range, the reviewer should determine if there are matrix effects, i.e. LCS data are within the QC limits but MS/MSD data are outside QC limit.

1. MS/MSD Recoveries and Precision Criteria

Data for MS and MSDs will not be present unless requested by the Region.
 Notify the Contract Laboratory Program Project Officer (CLP PO) if a field blank was used for the MS and MSD, unless designated as such by the Region.

NOTE: For a Matrix Spike that does not meet criteria, apply the action to only the field sample used to prepare the Matrix Spike sample. If it is clearly stated in the data validation materials that the samples were taken through incremental sampling or some other method guaranteeing the homogeneity of the sample group, then the entire sample group may be qualified.

List the %Rs, RPD of the compounds which do not meet the criteria.

Sample ID: JC19409-1 Matrix/Level: Aqueous

MS OR MSD	COMPOUND	% R	RPD	QC LIMITS	ACTION
<u> MS/MSD % recoveries and RPD within laboratory control limits except for the following: </u>					
<u> MS/MSD </u>	<u> Endrin aldehyde </u>	<u> 21/30 % </u>	<u> 34 </u>	<u> - 149/36 </u>	<u> No action </u>

Note: Analyte not detected in the samples analyzed, no qualification needed. MS/MSD results apply only to unspike sample.

Sample ID: JC19505-1 Matrix/Level: Soil

MS OR MSD	COMPOUND	% R	RPD	QC LIMITS	ACTION
<u> MS/MSD % recoveries and RPD within laboratory control limits except for the following: </u>					
<u> MS </u>	<u> 4,4'-DDE </u>	<u> 198 % </u>	<u> 124 </u>	<u> 10 - 155/49 </u>	<u> No action </u>
<u> MS </u>	<u> Endosulfan I </u>	<u> 178 % </u>	<u> 61 </u>	<u> 12 - 142/47 </u>	<u> No action </u>

Note: Analyte not detected in the samples analyzed, no qualification needed. MS/MSD results apply only to unspike sample.

DATA REVIEW WORKSHEETS

Action

No qualification of the data is necessary on MS and MSD data alone. However, using professional judgment, the validator may use the MS and MSD results in conjunction with other QC criteria and determine the need for some qualification of the data.

A separate worksheet should be used for each MS/MSD pair.

All criteria were met ☒
 Criteria were not met
 and/or see below _____

LABORATORY CONTROL SAMPLE (LCS) ANALYSIS

This data is generated to determine accuracy of the analytical method for various matrices.

1. LCS Recoveries Criteria

LCS Spike Compound	Recovery Limits (%)
gamma-BHC	50 – 120
Heptachlor epoxide	50 – 150
Dieldrin	30 – 130
4,4'-DDE	50 – 150
Endrin	50 – 120
Endosulfan sulfate	50 – 120
trans-Chlordane	30 – 130
Tetrachloro-m-xylene (surrogate)	30 – 150
Decachlorobiphenyl (surrogate)	30 – 150

LCS concentrations: 0.167 ug/L; 16.7 ug/kg

List the %R of compounds which do not meet the criteria

LCS ID	COMPOUND	% R	QC LIMIT

Action

The following guidance is suggested for qualifying sample data for which the associated LCS does not meet the required criteria.

- If the LCS recovery exceeds the upper acceptance limit, qualify detected target compounds as estimated (J). Do not qualify non-detected target compounds.
- If the LCS recovery is less than the lower acceptance limit, qualify detected target compounds as estimated (J) and non-detects as unusable (R).
- Use professional judgment to qualify data for compounds other than those compounds that are included in the LCS.
- Use professional judgment to qualify non-LCS compounds. Take into account the compound class, compound recovery efficiency, analytical problems associated with each compound, and comparability in the performance of the LCS compound to the non-LCS compound.
- If the LCS recovery is within allowable limits, no qualification of the data is necessary.

DATA REVIEW WORKSHEETS

2. Frequency Criteria:

Where LCS analyzed at the required frequency and for each matrix? Yes or No.

If no, the data may be affected. Use professional judgment to determine the severity of the effect and qualify data accordingly. Discuss any actions below and list the samples affected.

Note: Blank spike/blank spike duplicate analyzed for solid and aqueous matrices. % recoveries and RPD within laboratory control limits.

DATA REVIEW WORKSHEETS

All criteria were met _____
Criteria were not met _____
and/or see below N/A

FLORISIL CARTRIDGE PERFORMANCE CHECK

NOTE: Florisil cartridge cleanup is mandatory for all extracts.

Criteria

Is the Florisil cartridge performance check conducted at least once on each lot of cartridges used for sample cleanup or every 6 months, whichever is most frequent? Yes? or No? **N/A**

Criteria

Are the results for the Florisil Cartridge Performance Check solution included with the data package? Yes? or No? **N/A**

Note: If % criteria are not met, examine the raw data for the presence of polar interferences and use professional judgment in qualifying the data as follows:

Action:

- a. If the Percent Recovery is greater than 120% for any of the pesticide target compounds in the Florisil Cartridge Performance Check, qualify detected compounds as estimated (J). Do not qualify non-detected target compounds.
- b. If the Percent Recovery is greater than or equal to 80% and less than or equal to 120% for all the pesticide target compounds, no qualification of the data is necessary.
- c. If the Percent Recovery is greater than or equal to 10% and less than 80% for any of the pesticide target compounds in the Florisil Cartridge Performance Check, qualify detected target compounds as estimated (J) and non-detected target compounds as approximated (UJ).
- d. If the Percent Recovery is less than 10% for any of the pesticide target compounds in the Florisil Cartridge Performance Check, qualify detected compounds as estimated (J) and qualify non-detected target compounds as unusable (R).
- e. If the Percent Recovery of 2,4,5-trichlorophenol in the Florisil Cartridge Performance Check is greater than or equal to 5%, use professional judgment to qualify detected and non-detected target compounds, considering interference on the sample chromatogram.

Note: State in the Data Review Narrative potential effects on the sample data resulting from the Florisil Cartridge Performance Check analysis not yielding acceptable results.

Note: No information for florisil cartridge performance check included in data package. Florisil cartridge not used for sample extraction/clean-up. No qualification of the data performed, professional judgment.

All criteria were met N/A
Criteria were not met
and/or see below

GEL PERMEATION CHROMATOGRAPHY (GPC) PERFORMANCE CHECK

NOTE: GPC cleanup is mandatory for all soil samples.

If GPC criteria are not met, examine the raw data for the presence of high molecular weight contaminants; examine subsequent sample data for unusual peaks; and use professional judgment in qualifying the data. Notify the Contract Laboratory Program Project Officer (CLP PO) if the laboratory chooses to analyze samples under unacceptable GPC criteria.

Action:

- a. If the Percent Recovery is less than 10% for the pesticide compounds and surrogates during the GPC calibration check, the non-detected target compounds may be suspect, qualify detected compounds as estimated (J).
- b. If the Percent Recovery is less than 10% for the pesticide compounds and surrogates during the GPC calibration check, qualify all non-detected target compounds as unusable (R).
- c. If the Percent Recovery is greater than or equal to 10% and is less than 80% for any of the pesticide target compounds in the GPC calibration, qualify detected target compounds as estimated (J) and non-detected target compounds as approximated (UJ).
- d. If the Percent Recovery is greater than or equal to 80% and less than or equal to 120% for all the pesticide target compounds, no qualification of the data is necessary.
- e. If high recoveries (i.e., greater than 120%) were obtained for the pesticides and surrogates during the GPC calibration check, qualify detected compounds as estimated (J). Do not qualify non-detected target compounds.

Note: State in the Data Review Narrative potential effects on the sample data resulting from the GPC cleanup analyses not yielding acceptable results.

Note: **No information for performance of GPC cleanup included in data package. No qualification of the data performed, professional judgment.**

DATA REVIEW WORKSHEETS

All criteria were met X
Criteria were not met
and/or see below

TARGET COMPOUND IDENTIFICATION

Criteria:

1. Is Retention Times (RTs) of both of the surrogates and reported target compounds in each sample within the calculated RT Windows on both columns? Yes? or No?
2. Is the Tetrachloro-m-xylene (TCX) RT ± 0.05 minutes of the Mean RT (RT) determined from the initial calibration and Decachlorobiphenyl (DCB) within ± 0.10 minutes of the RT determined from the initial calibration? Yes? or No?
3. Is the Percent Difference (%D) for the detected mean concentrations of a pesticide target compound between the two Gas Chromatograph (GC) columns within the inclusive range of ± 25.0 %? Yes? or No?
4. When no analytes are identified in a sample; are the chromatograms from the analyses of the sample extract and the low-point standard of the initial calibration associated with those analyses on the same scaling factor? Yes? or No?
5. Does the chromatograms display the Single Component Pesticides (SCPs) detected in the sample and the largest peak of any multi-component analyte detected in the sample at less than full scale. Yes? or No?
6. If an extract is diluted; does the chromatogram display SCPs peaks between 10-100% of full scale, and multi-component analytes between 25-100% of full scale? Yes? or No? N/A
7. For any sample; does the baseline of the chromatogram return to below 50% of full scale before the elution time of alpha-BHC, and also return to below 25% of full scale after the elution time of alpha-BHC and before the elution time of DCB? Yes? or No?
8. If a chromatogram is replotted electronically to meet these requirements; is the scaling factor used displayed on the chromatogram, and both the initial chromatogram and the replotted chromatogram submitted in the data package. Yes? or No?

Action:

- a. If the qualitative criteria for both columns were not met, all target compounds that are reported as detected should be considered non-detected.
- b. Use professional judgment to assign an appropriate quantitation limit using the following guidance:
 - i. If the detected target compound peak was sufficiently outside the pesticide RT Window, the reported values may be a false positive and should be replaced with the sample Contract Required Quantitation Limits (CRQL) value.

DATA REVIEW WORKSHEETS

- ii. If the detected target compound peak poses an interference with potential detection of another target peak, the reported value should be considered and qualified as unusable (R).
- c. If the data reviewer identifies a peak in both GC column analyses that falls within the appropriate RT Windows, but was reported as a non-detect, the compound may be a false negative. Use professional judgment to decide if the compound should be included.

Note: State in the Data Review Narrative all conclusions made regarding target compound identification.

- d. If the Toxaphene peak RT windows determined from the calibration overlap with SCPs or chromatographic interferences, use professional judgment to qualify the data.
- e. If target compounds were detected on both GC columns, and the Percent Difference between the two results is greater than 25.0%, consider the potential for coelution and use professional judgment to decide whether a much larger concentration obtained on one column versus the other indicates the presence of an interfering compound. If an interfering compound is indicated, use professional judgment to determine how best to report, and if necessary, qualify the data according to these guidelines.
- f. If Toxaphene exhibits a marginal pattern-matching quality, use professional judgment to establish whether the differences are due to environmental "weathering" (i.e., degradation of the earlier eluting peaks relative to the later eluting peaks). If the presence of Toxaphene is strongly suggested, report results as presumptively present (N).

GAS CHROMATOGRAPH/MASS SPECTROMETER (GC/MS) CONFIRMATION

NOTE: This confirmation is not usually provided by the laboratory. In cases where it is provided, use professional judgment to determine if data qualified with "C" can be salvaged if it was previously qualified as unusable (R).

Action:

- a. If the quantitative criteria for both columns were met (≥ 5.0 ng/ μ L for SCPs and ≥ 125 ng/ μ L for Toxaphene), determine whether GC/MS confirmation was performed. If it was performed, qualify the data using the following guidance:
 - i. If GC/MS confirmation was not required because the quantitative criteria for both columns was not met, but it was still performed, use professional judgment when evaluating the data to decide whether the detect should be qualified with "C".
 - ii. If GC/MS confirmation was performed, but unsuccessful for a target compound detected by GC/ECD analysis, qualify those detects as "X".

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

COMPOUND QUANTITATION AND REPORTED CONTRACT REQUIRED QUANTITATION LIMITS (CRQLS)

The sample quantitation evaluation is to verify laboratory quantitation results. In the space below, please show a minimum of one sample calculation:

JC19473-1 Decachlorobiphenyl RF = 0.797

$$\begin{aligned} [] &= (53846125)(50)/(155.9 \times 10^6)(0.797) \\ &= 21.7 \text{ ppb} \quad \text{Ok} \end{aligned}$$

Action:

- If sample quantitation is different from the reported value, qualify result as unusable (R).
- When a sample is analyzed at more than one dilution, the lowest CRQLs are used unless a QC exceedance dictates the use of the higher CRQLs from the diluted sample.
- Replace concentrations that exceed the calibration range in the original analysis by crossing out the "E" and its corresponding value on the original reporting form and substituting the data from the diluted sample.
- Results between the MDL and CRQL should be qualified as estimated (J).
- Results less than the MDL should be reported at the CRQL and qualified (U). MDLs themselves are not reported.
- For non-aqueous samples, if the percent moisture is less than 70.0%, no qualification of the data is necessary. If the percent moisture is greater than or equal to 70.0% and less than 90.0%, qualify detects as estimated (J) and non-detects as approximated (UJ). If the percent moisture is greater than or equal to 90.0%, qualify detects as estimated (J) and non-detects as unusable (R) (see Table).

Percent Moisture Actions for Pesticide Analysis for Non-Aqueous Samples

Criteria	Action	
	Detected Associated Compounds	Non-detected Associated Compounds
% Moisture < 70.0	No qualification	
70.0 < % Moisture < 90.0	J	UJ
% Moisture > 90.0	J	R

DATA REVIEW WORKSHEETS

List samples which have $\leq 50\%$ solids

Note: If any discrepancies are found, the Region's designated representative may contact the laboratory to obtain additional information that could resolve any differences. If a discrepancy remains unresolved, the reviewer must use professional judgment to decide which value is the most accurate. Under these circumstances, the reviewer may determine that qualification of data is warranted. Note in the Data Review Narrative a description of the reasons for data qualification and the qualification that is applied to the data.

Dilution performed

SAMPLE ID	DILUTION FACTOR	REASON FOR DILUTION

DATA REVIEW WORKSHEETS

All criteria were met NA
 Criteria were not met
 and/or see below _____

FIELD DUPLICATE PRECISION

NOTE: In the absence of QAPP guidance for validating data from field duplicates, the following action will be taken.

Field duplicate samples may be taken and analyzed as an indication of overall precision. These analyses measure both field and lab precision; therefore, the results may have more variability than laboratory duplicates which only laboratory performance. It is also expected that soil duplicate results will have a greater variance than water matrices due to difficulties associated with collecting identical field duplicate samples. Identify which samples within the data package are field duplicates. Estimate the relative percent difference (RPD) between the values for each compound. If large RPDs (> 50%) is observed, confirm identification of samples and note difference in the executive summary.

Sample IDs: _____ - _____

Matrix: _____ - _____

COMPOUND	SQL ug/L	SAMPLE CONC.	DUPLICATE CONC.	RPD	ACTION
No field/laboratory duplicate analyzed with this data package. MS/MSD or LCS/LCSD % recoveries RPD used to assess precision. RPD within the required criteria of < 50 %.					

Actions:

a. Qualify as estimated positive results (J) and nondetects (UJ) for the compound that exceeded the above criteria. For organics, only the sample and duplicate will be qualified.

b. If an RPD cannot be calculated because one or both of the sample results is not detected, the following actions apply:

- i. If one sample result is not detected and the other is greater than 5x the SQL qualify (J/UJ).
- ii. If one sample value is not detected and the other is greater than 5x the SQL and the SQLs for the sample and duplicate are significantly different, use professional judgment to determine if qualification is appropriate.
- iii. If one sample value is not detected and the other is less than 5x, use professional judgment to determine if qualification is appropriate.
- iv. If both sample and duplicate results are not detected, no action is needed.

OVERALL ASSESSMENT OF DATA

Action:

1. Use professional judgment to determine if there is any need to qualify data which were not qualified based on the Quality Control (QC) criteria previously discussed.
2. Write a brief narrative to give the user an indication of the analytical limitations of the data.

Note: The Contract Laboratory Program Project Officer (CLP PO) must be informed if any inconsistency of the data with the Sample Delivery Group (SDG) Narrative. If sufficient information on the intended use and required quality of the data is available, the reviewer should include their assessment of the usability of the data within the given context. This may be used as part of a formal Data Quality Assessment (DQA).

Overall assessment of the data: Results are valid; the data can be used for decision making purposes.