

CETIFICATION

SDG No: JC19644 Laboratory: Accutest, New Jersey
 Site: BMS, Building 5 Area, PR Matrix: Soil/Groundwater
 Humacao, PR

SUMMARY: Soil and groundwater samples (Table 1) were collected on the BMSMC facility – Building 5 Area. The BMSMC facility is located in Humacao, PR. Samples were taken May 3-4, 2016 and were analyzed in Accutest Laboratory of Dayton, New Jersey for the ABN TCL Special List (1,4-Dioxane and Naphthalene were analyzed following the SIM technique); TCL pesticides list; and for low molecular weight alcohols (LMWA) the results were reported under SDG No.: JC19644. Results were validated using the latest validation guidelines (July, 2015) of the EPA Hazardous Waste Support Section. The analyses performed are shown in Table 1. Individual data review worksheets are enclosed for each target analyte group. The data sample organic data samples summary form shows for analytes results that were qualified.

In summary the results are valid and can be used for decision taking purposes.

Table 1. Samples analyzed and analysis performed

SAMPLE ID	SAMPLE DESCRIPTION	MATRIX	ANALYSIS PERFORMED
JC19644-1	RA-16 GWD	Groundwater	ABN TCL special list; pesticides TCL list; LMWA
JC19644-2	S-42S	Groundwater	ABN TCL special list; pesticides TCL list; LMWA
JC19644-3	RA-11 GWS	Groundwater	ABN TCL special list; pesticides TCL list; LMWA
JC19644-4	RA-11 (10-11)	Soil	ABN TCL special list; pesticides TCL list; LMWA
JC19644-5	RA-10 (5.5-6.5)	Soil	ABN TCL special list; pesticides TCL list; LMWA

Reviewer Name: Rafael Infante
 Chemist License 1888

Signature:



Date:

May 30, 2016



SGS Accutest

Report of Analysis

Page 1 of 3

Client Sample ID:	RA-16 GWD	Date Sampled:	05/03/16
Lab Sample ID:	JC19644-1	Date Received:	05/05/16
Matrix:	AQ - Ground Water	Percent Solids:	n/a
Method:	SW846 8270D SW846 3510C		
Project:	BMSMC, Building 5 Area, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	F156860.D	1	05/06/16	LK	05/05/16	OP93662	EF6604
Run #2							

Run #	Initial Volume	Final Volume
Run #1	900 ml	1.0 ml
Run #2		

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units	Q
95-57-8	2-Chlorophenol	ND	5.6	0.91	ug/l	
59-50-7	4-Chloro-3-methyl phenol	ND	5.6	0.99	ug/l	
120-83-2	2,4-Dichlorophenol	ND	2.2	1.4	ug/l	
105-67-9	2,4-Dimethylphenol	ND	5.6	2.7	ug/l	
51-28-5	2,4-Dinitrophenol	ND	11	1.7	ug/l	
534-52-1	4,6-Dinitro-o-cresol	ND	5.6	1.4	ug/l	
95-48-7	2-Methylphenol	ND	2.2	0.99	ug/l	
	3&4-Methylphenol	ND	2.2	0.98	ug/l	
88-75-5	2-Nitrophenol	ND	5.6	1.1	ug/l	
100-02-7	4-Nitrophenol	ND	11	1.3	ug/l	
87-86-5	Pentachlorophenol	ND	5.6	1.5	ug/l	
108-95-2	Phenol	ND	2.2	0.44	ug/l	
58-90-2	2,3,4,6-Tetrachlorophenol	ND	5.6	1.6	ug/l	
95-95-4	2,4,5-Trichlorophenol	ND	5.6	1.5	ug/l	
88-06-2	2,4,6-Trichlorophenol	ND	5.6	1.0	ug/l	
83-32-9	Acenaphthene	ND	1.1	0.21	ug/l	
208-96-8	Acenaphthylene	ND	1.1	0.15	ug/l	
98-86-2	Acetophenone	ND	2.2	0.23	ug/l	
120-12-7	Anthracene	ND	1.1	0.23	ug/l	
1912-24-9	Atrazine	ND	2.2	0.50	ug/l	
100-52-7	Benzaldehyde	ND	5.6	0.32	ug/l	
56-55-3	Benzo(a)anthracene	ND	1.1	0.23	ug/l	
50-32-8	Benzo(a)pyrene	ND	1.1	0.24	ug/l	
205-99-2	Benzo(b)fluoranthene	ND	1.1	0.23	ug/l	
191-24-2	Benzo(g,h,i)perylene	ND	1.1	0.38	ug/l	
207-08-9	Benzo(k)fluoranthene	ND	1.1	0.23	ug/l	
101-55-3	4-Bromophenyl phenyl ether	ND	2.2	0.45	ug/l	
85-68-7	Butyl benzyl phthalate	ND	2.2	0.51	ug/l	
92-52-4	1,1'-Biphenyl	ND	1.1	0.24	ug/l	
91-58-7	2-Chloronaphthalene	ND	2.2	0.26	ug/l	
106-47-8	4-Chloroaniline	ND	5.6	0.38	ug/l	
86-74-8	Carbazole	ND	1.1	0.25	ug/l	

ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound



Report of Analysis

Client Sample ID: RA-16 GWD
 Lab Sample ID: JC19644-1
 Matrix: AQ - Ground Water
 Method: SW846 8270D SW846 3510C
 Project: BMSMC, Building 5 Area, PR

Date Sampled: 05/03/16
 Date Received: 05/05/16
 Percent Solids: n/a

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units	Q
105-60-2	Caprolactam	ND	2.2	0.72	ug/l	
218-01-9	Chrysene	ND	1.1	0.20	ug/l	
111-91-1	bis(2-Chloroethoxy)methane	ND	2.2	0.31	ug/l	
111-44-4	bis(2-Chloroethyl)ether	ND	2.2	0.28	ug/l	
108-60-1	bis(2-Chloroisopropyl)ether	ND	2.2	0.45	ug/l	
7005-72-3	4-Chlorophenyl phenyl ether	ND	2.2	0.41	ug/l	
121-14-2	2,4-Dinitrotoluene	ND	1.1	0.61	ug/l	
606-20-2	2,6-Dinitrotoluene	ND	1.1	0.53	ug/l	
91-94-1	3,3'-Dichlorobenzidine	ND	2.2	0.56	ug/l	
123-91-1	1,4-Dioxane	26.7	1.1	0.73	ug/l	
53-70-3	Dibenzo(a,h)anthracene	ND	1.1	0.37	ug/l	
132-64-9	Dibenzofuran	ND	5.6	0.24	ug/l	
84-74-2	Di-n-butyl phthalate	ND	2.2	0.55	ug/l	
117-84-0	Di-n-octyl phthalate	ND	2.2	0.26	ug/l	
84-66-2	Diethyl phthalate	ND	2.2	0.29	ug/l	
131-11-3	Dimethyl phthalate	ND	2.2	0.24	ug/l	
117-81-7	bis(2-Ethylhexyl)phthalate	ND	2.2	1.8	ug/l	
206-44-0	Fluoranthene	ND	1.1	0.19	ug/l	
86-73-7	Fluorene	ND	1.1	0.19	ug/l	
118-74-1	Hexachlorobenzene	ND	1.1	0.36	ug/l	
87-68-3	Hexachlorobutadiene	ND	1.1	0.55	ug/l	
77-47-4	Hexachlorocyclopentadiene	ND	11	3.1	ug/l	
67-72-1	Hexachloroethane	ND	2.2	0.43	ug/l	
193-39-5	Indeno(1,2,3-cd)pyrene	ND	1.1	0.37	ug/l	
78-59-1	Isophorone	ND	2.2	0.31	ug/l	
90-12-0	1-Methylnaphthalene	ND	1.1	0.29	ug/l	
91-57-6	2-Methylnaphthalene	ND	1.1	0.23	ug/l	
88-74-1	2-Nitroaniline	ND	5.6	0.31	ug/l	
99-09-2	3-Nitroaniline	ND	5.6	0.43	ug/l	
100-01-6	4-Nitroaniline	ND	5.6	0.49	ug/l	
98-95-3	Nitrobenzene	ND	2.2	0.71	ug/l	
621-64-7	N-Nitroso-di-n-propylamine	ND	2.2	0.53	ug/l	
86-30-6	N-Nitrosodiphenylamine	ND	5.6	0.25	ug/l	
85-01-8	Phenanthrene	ND	1.1	0.19	ug/l	
129-00-0	Pyrene	ND	1.1	0.24	ug/l	
95-94-3	1,2,4,5-Tetrachlorobenzene	ND	2.2	0.41	ug/l	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
367-12-4	2-Fluorophenol	42%		14-88%

ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound



Report of Analysis

Client Sample ID: RA-16 GWD	
Lab Sample ID: JC19644-1	Date Sampled: 05/03/16
Matrix: AQ - Ground Water	Date Received: 05/05/16
Method: SW846 8270D SW846 3510C	Percent Solids: n/a
Project: BMSMC, Building 5 Area, PR	

ABN TCL Special List

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
4165-62-2	Phenol-d5	31%		10-110%
118-79-6	2,4,6-Tribromophenol	93%		39-149%
4165-60-0	Nitrobenzene-d5	78%		32-128%
321-60-8	2-Fluorobiphenyl	68%		35-119%
1718-51-0	Terphenyl-d14	73%		10-126%



ND = Not detected MDL = Method Detection Limit
RL = Reporting Limit
E = Indicates value exceeds calibration range

J = Indicates an estimated value
B = Indicates analyte found in associated method blank
N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID:	RA-16 GWD	Date Sampled:	05/03/16
Lab Sample ID:	JC19644-1	Date Received:	05/05/16
Matrix:	AQ - Ground Water	Percent Solids:	n/a
Method:	SW846 8270D BY SIM SW846 3510C		
Project:	BMSMC, Building 5 Area, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	3M61128.D	1	05/06/16	LK	05/05/16	OP93662A	E3M2874
Run #2							

Run #	Initial Volume	Final Volume
Run #1	900 ml	1.0 ml
Run #2		

CAS No.	Compound	Result	RL	MDL	Units	Q
91-20-3	Naphthalene	ND	0.11	0.033	ug/l	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
4165-60-0	Nitrobenzene-d5	70%		24-125%
321-60-8	2-Fluorobiphenyl	72%		19-127%
1718-51-0	Terphenyl-d14	84%		10-119%



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID: RA-16 GWD		Date Sampled: 05/03/16
Lab Sample ID: JC19644-1		Date Received: 05/05/16
Matrix: AQ - Ground Water		Percent Solids: n/a
Method: SW846-8015C (DAD)		
Project: BSMC, Building 5 Area, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	GH104963.D	1	05/16/16	XPL	n/a	n/a	GGH5285
Run #2							

Low Molecular Alcohol List

CAS No.	Compound	Result	RL	MDL	Units	Q
64-17-5	Ethanol	ND	100	55	ug/l	
78-83-1	Isobutyl Alcohol	ND	100	36	ug/l	
67-63-0	Isopropyl Alcohol	ND	100	68	ug/l	
71-23-8	n-Propyl Alcohol	ND	100	43	ug/l	
71-36-3	n-Butyl Alcohol	ND	100	87	ug/l	
78-92-2	sec-Butyl Alcohol	ND	100	66	ug/l	
67-56-1	Methanol	ND	200	71	ug/l	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
111-27-3	Hexanol	104%		56-145%



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID: RA-16 GWD
 Lab Sample ID: JC19644-1
 Matrix: AQ - Ground Water
 Method: SW846 8081B SW846 3510C
 Project: BMSMC, Building 5 Area, PR

Date Sampled: 05/03/16
 Date Received: 05/05/16
 Percent Solids: n/a

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	4G67960.D	1	05/06/16	BP	05/05/16	OP93665	G4G1782
Run #2							

Run #	Initial Volume	Final Volume
Run #1	300 ml	2.0 ml
Run #2		

Pesticide TCL List

CAS No.	Compound	Result	RL	MDL	Units	Q
309-00-2	Aldrin	ND	0.0067	0.0040	ug/l	
319-84-6	alpha-BHC	ND	0.0067	0.0040	ug/l	
319-85-7	beta-BHC	ND	0.0067	0.0038	ug/l	
319-86-8	delta-BHC	ND	0.0067	0.0030	ug/l	
58-89-9	gamma-BHC (Lindane)	ND	0.0067	0.0019	ug/l	
5103-71-9	alpha-Chlordane	ND	0.0067	0.0031	ug/l	
5103-74-2	gamma-Chlordane	ND	0.0067	0.0031	ug/l	
60-57-1	Dieldrin	ND	0.0067	0.0024	ug/l	
72-54-8	4,4'-DDD	ND	0.0067	0.0025	ug/l	
72-55-9	4,4'-DDE	ND	0.0067	0.0041	ug/l	
50-29-3	4,4'-DDT	ND	0.0067	0.0033	ug/l	
72-20-8	Endrin	ND	0.0067	0.0034	ug/l	
1031-07-8	Endosulfan sulfate	ND	0.0067	0.0035	ug/l	
7421-93-4	Endrin aldehyde	ND	0.0067	0.0034	ug/l	
53494-70-5	Endrin ketone	ND	0.0067	0.0034	ug/l	
959-98-8	Endosulfan-I	ND	0.0067	0.0033	ug/l	
33213-65-9	Endosulfan-II	ND	0.0067	0.0029	ug/l	
76-44-8	Heptachlor	ND	0.0067	0.0025	ug/l	
1024-57-3	Heptachlor epoxide	ND	0.0067	0.0044	ug/l	
72-43-5	Methoxychlor	ND	0.013	0.0038	ug/l	
8001-35-2	Toxaphene	ND	0.17	0.12	ug/l	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
877-09-8	Tetrachloro-m-xylene	67%		26-132%
877-09-8	Tetrachloro-m-xylene	68%		26-132%
2051-24-3	Decachlorobiphenyl	107%		10-118%
2051-24-3	Decachlorobiphenyl	95%		10-118%



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 3

Client Sample ID:	S-425	Date Sampled:	05/03/16
Lab Sample ID:	JC19644-2	Date Received:	05/05/16
Matrix:	AQ - Ground Water	Percent Solids:	n/a
Method:	SW846 8270D SW846 3510C		
Project:	BMSMC, Building 5 Area, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	F156861.D	1	05/06/16	LK	05/05/16	OP93662	EF6604
Run #2	F156880.D	50	05/06/16	LK	05/05/16	OP93662	EF6604

Run #	Initial Volume	Final Volume
Run #1	1000 ml	1.0 ml
Run #2	1000 ml	1.0 ml

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units	Q
95-57-8	2-Chlorophenol	ND	5.0	0.82	ug/l	
59-50-7	4-Chloro-3-methyl phenol	ND	5.0	0.89	ug/l	
120-83-2	2,4-Dichlorophenol	ND	2.0	1.3	ug/l	
105-67-9	2,4-Dimethylphenol	ND	5.0	2.4	ug/l	
51-28-5	2,4-Dinitrophenol	ND	10	1.6	ug/l	
534-52-1	4,6-Dinitro-o-cresol	ND	5.0	1.3	ug/l	
95-48-7	2-Methylphenol	ND	2.0	0.89	ug/l	
	3&4-Methylphenol	ND	2.0	0.88	ug/l	
88-75-5	2-Nitrophenol	ND	5.0	0.96	ug/l	
100-02-7	4-Nitrophenol	ND	10	1.2	ug/l	
87-86-5	Pentachlorophenol	ND	5.0	1.4	ug/l	
108-95-2	Phenol	ND	2.0	0.39	ug/l	
58-90-2	2,3,4,6-Tetrachlorophenol	ND	5.0	1.5	ug/l	
95-95-4	2,4,5-Trichlorophenol	ND	5.0	1.3	ug/l	
88-06-2	2,4,6-Trichlorophenol	ND	5.0	0.92	ug/l	
83-32-9	Acenaphthene	ND	1.0	0.19	ug/l	
208-96-8	Acenaphthylene	ND	1.0	0.14	ug/l	
98-86-2	Acetophenone	ND	2.0	0.21	ug/l	
120-12-7	Anthracene	ND	1.0	0.21	ug/l	
1912-24-9	Atrazine	ND	2.0	0.45	ug/l	
100-52-7	Benzaldehyde	ND	5.0	0.29	ug/l	
56-55-3	Benzo(a)anthracene	ND	1.0	0.20	ug/l	
50-32-8	Benzo(a)pyrene	ND	1.0	0.21	ug/l	
205-99-2	Benzo(b)fluoranthene	ND	1.0	0.21	ug/l	
191-24-2	Benzo(g,h,i)perylene	ND	1.0	0.34	ug/l	
207-08-9	Benzo(k)fluoranthene	ND	1.0	0.21	ug/l	
101-55-3	4-Bromophenyl phenyl ether	ND	2.0	0.40	ug/l	
85-68-7	Butyl benzyl phthalate	ND	2.0	0.46	ug/l	
92-52-4	1,1'-Biphenyl	ND	1.0	0.21	ug/l	
91-58-7	2-Chloronaphthalene	ND	2.0	0.24	ug/l	
106-47-8	4-Chloroaniline	ND	5.0	0.34	ug/l	
86-74-8	Carbazole	ND	1.0	0.23	ug/l	



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

Report of Analysis

Client Sample ID: S-42S
 Lab Sample ID: JC19644-2
 Matrix: AQ - Ground Water
 Method: SW846 8270D SW846 3510C
 Project: BSMC, Building 5 Area, PR

Date Sampled: 05/03/16
 Date Received: 05/05/16
 Percent Solids: n/a

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units	Q
105-60-2	Caprolactam	ND	2.0	0.65	ug/l	
218-01-9	Chrysene	ND	1.0	0.18	ug/l	
111-91-1	bis(2-Chloroethoxy)methane	ND	2.0	0.28	ug/l	
111-44-4	bis(2-Chloroethyl)ether	ND	2.0	0.25	ug/l	
108-60-1	bis(2-Chloroisopropyl)ether	ND	2.0	0.40	ug/l	
7005-72-3	4-Chlorophenyl phenyl ether	ND	2.0	0.37	ug/l	
121-14-2	2,4-Dinitrotoluene	ND	1.0	0.55	ug/l	
606-20-2	2,6-Dinitrotoluene	ND	1.0	0.48	ug/l	
91-94-1	3,3'-Dichlorobenzidine	ND	2.0	0.51	ug/l	
123-91-1	1,4-Dioxane	3890 ^a	50	33	ug/l	
53-70-3	Dibenzo(a,h)anthracene	ND	1.0	0.33	ug/l	
132-64-9	Dibenzofuran	ND	5.0	0.22	ug/l	
84-74-2	Di-n-butyl phthalate	ND	2.0	0.50	ug/l	
117-84-0	Di-n-octyl phthalate	ND	2.0	0.23	ug/l	
84-66-2	Diethyl phthalate	ND	2.0	0.26	ug/l	
131-11-3	Dimethyl phthalate	ND	2.0	0.22	ug/l	
117-81-7	bis(2-Ethylhexyl)phthalate	ND	2.0	1.7	ug/l	
206-44-0	Fluoranthene	ND	1.0	0.17	ug/l	
86-73-7	Fluorene	ND	1.0	0.17	ug/l	
118-74-1	Hexachlorobenzene	ND	1.0	0.33	ug/l	
87-68-3	Hexachlorobutadiene	ND	1.0	0.49	ug/l	
77-47-4	Hexachlorocyclopentadiene	ND	10	2.8	ug/l	
67-72-1	Hexachloroethane	ND	2.0	0.39	ug/l	
193-39-5	Indeno(1,2,3-cd)pyrene	ND	1.0	0.33	ug/l	
78-59-1	Isophorone	ND	2.0	0.28	ug/l	
90-12-0	1-Methylnaphthalene	ND	1.0	0.26	ug/l	
91-57-6	2-Methylnaphthalene	ND	1.0	0.21	ug/l	
88-74-4	2-Nitroaniline	ND	5.0	0.28	ug/l	
99-09-2	3-Nitroaniline	ND	5.0	0.39	ug/l	
100-01-6	4-Nitroaniline	ND	5.0	0.44	ug/l	
98-95-3	Nitrobenzene	ND	2.0	0.64	ug/l	
621-64-7	N-Nitroso-di-n-propylamine	ND	2.0	0.48	ug/l	
86-30-6	N-Nitrosodiphenylamine	ND	5.0	0.22	ug/l	
85-01-8	Phenanthrene	ND	1.0	0.18	ug/l	
129-00-0	Pyrene	ND	1.0	0.22	ug/l	
95-94-3	1,2,4,5-Tetrachlorobenzene	ND	2.0	0.37	ug/l	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
367-12-4	2-Fluorophenol	41%	0% ^b	14-88%

ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound



Report of Analysis

Client Sample ID:	S-42S	Date Sampled:	05/03/16
Lab Sample ID:	JC19644-2	Date Received:	05/05/16
Matrix:	AQ - Ground Water	Percent Solids:	n/a
Method:	SW846 8270D SW846 3510C		
Project:	BMSMC, Building 5 Area, PR		

ABN TCL Special List

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
4165-62-2	Phenol-d5	27%	0% ^b	10-110%
118-79-6	2,4,6-Tribromophenol	98%	0% ^b	39-149%
4165-60-0	Nitrobenzene-d5	74%	70%	32-128%
321-60-8	2-Fluorobiphenyl	67%	63%	35-119%
1718-51-0	Terphenyl-d14	82%	71%	10-126%

(a) Result is from Run# 2

(b) Outside control limits due to dilution.



ND = Not detected MDL = Method Detection Limit
RL = Reporting Limit
E = Indicates value exceeds calibration range

J = Indicates an estimated value
B = Indicates analyte found in associated method blank
N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID:	S-42S	Date Sampled:	05/03/16
Lab Sample ID:	JC19644-2	Date Received:	05/05/16
Matrix:	AQ - Ground Water	Percent Solids:	n/a
Method:	SW846 8270D BY SIM SW846 3510C		
Project:	BMSMC, Building 5 Area, PR		

	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	3M61129.D	1	05/06/16	LK	05/05/16	OP93662A	E3M2874
Run #2							

	Initial Volume	Final Volume
Run #1	1000 ml	1.0 ml
Run #2		

CAS No.	Compound	Result	RL	MDL	Units	Q
91-20-3	Naphthalene	ND	0.10	0.029	ug/l	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
4165-60-0	Nitrobenzene-d5	66%		24-125%
321-60-8	2-Fluorobiphenyl	67%		19-127%
1718-51-0	Terphenyl-d14	91%		10-119%



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID:	S-42S	Date Sampled:	05/03/16
Lab Sample ID:	JC19644-2	Date Received:	05/05/16
Matrix:	AQ - Ground Water	Percent Solids:	n/a
Method:	SW846-8015C (DAD)		
Project:	BMSMC, Building 5 Area, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	GH104966.D	1	05/16/16	XPL	n/a	n/a	GGH5285
Run #2							

Low Molecular Alcohol List

CAS No.	Compound	Result	RL	MDL	Units	Q
64-17-5	Ethanol	ND	100	55	ug/l	
78-83-1	Isobutyl Alcohol	ND	100	36	ug/l	
67-63-0	Isopropyl Alcohol	ND	100	68	ug/l	
71-23-8	n-Propyl Alcohol	ND	100	43	ug/l	
71-36-3	n-Butyl Alcohol	ND	100	87	ug/l	
78-92-2	sec-Butyl Alcohol	ND	100	66	ug/l	
67-56-1	Methanol	ND	200	71	ug/l	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
111-27-3	Hexanol	90%		56-145%



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID:	S-42S	Date Sampled:	05/03/16
Lab Sample ID:	JC19644-2	Date Received:	05/05/16
Matrix:	AQ - Ground Water	Percent Solids:	n/a
Method:	SW846 8081B SW846 3510C		
Project:	BMSMC, Building 5 Area, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	4G67961.D	1	05/06/16	BP	05/05/16	OP93665	G4G1782
Run #2							

Run #	Initial Volume	Final Volume
Run #1	300 ml	2.0 ml
Run #2		

Pesticide TCL List

CAS No.	Compound	Result	RL	MDL	Units	Q
309-00-2	Aldrin	ND	0.0067	0.0040	ug/l	
319-84-6	alpha-BHC	ND	0.0067	0.0040	ug/l	
319-85-7	beta-BHC	ND	0.0067	0.0038	ug/l	
319-86-8	delta-BHC	ND	0.0067	0.0030	ug/l	
58-89-9	gamma-BHC (Lindane)	ND	0.0067	0.0019	ug/l	
5103-71-9	alpha-Chlordane	ND	0.0067	0.0031	ug/l	
5103-74-2	gamma-Chlordane	ND	0.0067	0.0031	ug/l	
60-57-1	Dieldrin	ND	0.0067	0.0024	ug/l	
72-54-8	4,4'-DDD	ND	0.0067	0.0025	ug/l	
72-55-9	4,4'-DDE	ND	0.0067	0.0041	ug/l	
50-29-3	4,4'-DDT	ND	0.0067	0.0033	ug/l	
72-20-8	Endrin	ND	0.0067	0.0034	ug/l	
1031-07-8	Endosulfan sulfate	ND	0.0067	0.0035	ug/l	
7421-93-4	Endrin aldehyde	ND	0.0067	0.0034	ug/l	
53494-70-5	Endrin ketone	ND	0.0067	0.0034	ug/l	
959-98-8	Endosulfan-I	ND	0.0067	0.0033	ug/l	
33213-65-9	Endosulfan-II	ND	0.0067	0.0029	ug/l	
76-44-8	Heptachlor	ND	0.0067	0.0025	ug/l	
1024-57-3	Heptachlor epoxide	ND	0.0067	0.0044	ug/l	
72-43-5	Methoxychlor	ND	0.013	0.0038	ug/l	
8001-35-2	Toxaphene	ND	0.17	0.12	ug/l	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
877-09-8	Tetrachloro-m-xylene	62%		26-132%
877-09-8	Tetrachloro-m-xylene	57%		26-132%
2051-24-3	Decachlorobiphenyl	107%		10-118%
2051-24-3	Decachlorobiphenyl	88%		10-118%



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 3

Client Sample ID: RA-11 GWS
 Lab Sample ID: JC19644-3
 Matrix: AQ - Ground Water
 Method: SW846 8270D SW846 3510C
 Project: BSMC, Building 5 Area, PR

Date Sampled: 05/03/16
 Date Received: 05/05/16
 Percent Solids: n/a

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	F156862.D	1	05/06/16	LK	05/05/16	OP93662	EF6604
Run #2							

Run #	Initial Volume	Final Volume
Run #1	900 ml	1.0 ml
Run #2		

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units	Q
95-57-8	2-Chlorophenol	ND	5.6	0.91	ug/l	
59-50-7	4-Chloro-3-methyl phenol	ND	5.6	0.99	ug/l	
120-83-2	2,4-Dichlorophenol	ND	2.2	1.4	ug/l	
105-67-9	2,4-Dimethylphenol	ND	5.6	2.7	ug/l	
51-28-5	2,4-Dinitrophenol	ND	11	1.7	ug/l	
534-52-1	4,6-Dinitro-o-cresol	ND	5.6	1.4	ug/l	
95-48-7	2-Methylphenol	ND	2.2	0.99	ug/l	
	3&4-Methylphenol	ND	2.2	0.98	ug/l	
88-75-5	2-Nitrophenol	ND	5.6	1.1	ug/l	
100-02-7	4-Nitrophenol	ND	11	1.3	ug/l	
87-86-5	Pentachlorophenol	ND	5.6	1.5	ug/l	
108-95-2	Phenol	ND	2.2	0.44	ug/l	
58-90-2	2,3,4,6-Tetrachlorophenol	ND	5.6	1.6	ug/l	
95-95-4	2,4,5-Trichlorophenol	ND	5.6	1.5	ug/l	
88-06-2	2,4,6-Trichlorophenol	ND	5.6	1.0	ug/l	
83-32-9	Acenaphthene	ND	1.1	0.21	ug/l	
208-96-8	Acenaphthylene	ND	1.1	0.15	ug/l	
98-86-2	Acetophenone	ND	2.2	0.23	ug/l	
120-12-7	Anthracene	ND	1.1	0.23	ug/l	
1912-24-9	Atrazine	ND	2.2	0.50	ug/l	
100-52-7	Benzaldehyde	ND	5.6	0.32	ug/l	
56-55-3	Benzo(a)anthracene	ND	1.1	0.23	ug/l	
50-32-8	Benzo(a)pyrene	ND	1.1	0.24	ug/l	
205-99-2	Benzo(b)fluoranthene	ND	1.1	0.23	ug/l	
191-24-2	Benzo(g,h,i)perylene	ND	1.1	0.38	ug/l	
207-08-9	Benzo(k)fluoranthene	ND	1.1	0.23	ug/l	
101-55-3	4-Bromophenyl phenyl ether	ND	2.2	0.45	ug/l	
85-68-7	Butyl benzyl phthalate	ND	2.2	0.51	ug/l	
92-52-4	1,1'-Biphenyl	ND	1.1	0.24	ug/l	
91-58-7	2-Chloronaphthalene	ND	2.2	0.26	ug/l	
106-47-8	4-Chloroaniline	ND	5.6	0.38	ug/l	
86-74-8	Carbazole	ND	1.1	0.25	ug/l	



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

Report of Analysis

Client Sample ID: RA-11 GWS
 Lab Sample ID: JC19644-3
 Matrix: AQ - Ground Water
 Method: SW846 8270D SW846 3510C
 Project: BMSMC, Building 5 Area, PR

Date Sampled: 05/03/16
 Date Received: 05/05/16
 Percent Solids: n/a

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units	Q
105-60-2	Caprolactam	ND	2.2	0.72	ug/l	
218-01-9	Chrysene	ND	1.1	0.20	ug/l	
111-91-1	bis(2-Chloroethoxy)methane	ND	2.2	0.31	ug/l	
111-44-4	bis(2-Chloroethyl)ether	ND	2.2	0.28	ug/l	
108-60-1	bis(2-Chloroisopropyl)ether	ND	2.2	0.45	ug/l	
7005-72-3	4-Chlorophenyl phenyl ether	ND	2.2	0.41	ug/l	
121-14-2	2,4-Dinitrotoluene	ND	1.1	0.61	ug/l	
606-20-2	2,6-Dinitrotoluene	ND	1.1	0.53	ug/l	
91-94-1	3,3'-Dichlorobenzidine	ND	2.2	0.56	ug/l	
123-91-1	1,4-Dioxane	67.1	1.1	0.73	ug/l	
53-70-3	Dibenzo(a,h)anthracene	ND	1.1	0.37	ug/l	
132-64-9	Dibenzofuran	ND	5.6	0.24	ug/l	
84-74-2	Di-n-butyl phthalate	ND	2.2	0.55	ug/l	
117-84-0	Di-n-octyl phthalate	ND	2.2	0.26	ug/l	
84-66-2	Diethyl phthalate	ND	2.2	0.29	ug/l	
131-11-3	Dimethyl phthalate	ND	2.2	0.24	ug/l	
117-81-7	bis(2-Ethylhexyl)phthalate	ND	2.2	1.8	ug/l	
206-44-0	Fluoranthene	ND	1.1	0.19	ug/l	
86-73-7	Fluorene	ND	1.1	0.19	ug/l	
118-74-1	Hexachlorobenzene	ND	1.1	0.36	ug/l	
87-68-3	Hexachlorobutadiene	ND	1.1	0.55	ug/l	
77-47-4	Hexachlorocyclopentadiene	ND	11	3.1	ug/l	
67-72-1	Hexachloroethane	ND	2.2	0.43	ug/l	
193-39-5	Indeno(1,2,3-cd)pyrene	ND	1.1	0.37	ug/l	
78-59-1	Isophorone	ND	2.2	0.31	ug/l	
90-12-0	1-Methylnaphthalene	ND	1.1	0.29	ug/l	
91-57-6	2-Methylnaphthalene	ND	1.1	0.23	ug/l	
88-74-1	2-Nitroaniline	ND	5.6	0.31	ug/l	
99-09-2	3-Nitroaniline	ND	5.6	0.43	ug/l	
100-01-6	4-Nitroaniline	ND	5.6	0.49	ug/l	
98-95-3	Nitrobenzene	ND	2.2	0.71	ug/l	
621-64-7	N-Nitroso-di-n-propylamine	ND	2.2	0.53	ug/l	
86-30-6	N-Nitrosodiphenylamine	ND	5.6	0.25	ug/l	
85-01-8	Phenanthrene	ND	1.1	0.19	ug/l	
129-00-0	Pyrene	ND	1.1	0.24	ug/l	
95-94-3	1,2,4,5-Tetrachlorobenzene	ND	2.2	0.41	ug/l	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
367-12-4	2-Fluorophenol	42%		14-88%

ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound



Report of Analysis

Client Sample ID: RA-11 GWS
Lab Sample ID: JC19644-3
Matrix: AQ - Ground Water
Method: SW846 8270D SW846 3510C
Project: BSMC, Building 5 Area, PR

Date Sampled: 05/03/16
Date Received: 05/05/16
Percent Solids: n/a

ABN TCL Special List

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
4165-62-2	Phenol-d5	29%		10-110%
118-79-6	2,4,6-Tribromophenol	89%		39-149%
4165-60-0	Nitrobenzene-d5	73%		32-128%
321-60-8	2-Fluorobiphenyl	64%		35-119%
1718-51-0	Terphenyl-d14	64%		10-126%



ND = Not detected MDL = Method Detection Limit
RL = Reporting Limit
E = Indicates value exceeds calibration range

J = Indicates an estimated value
B = Indicates analyte found in associated method blank
N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID:	RA-11 GWS	Date Sampled:	05/03/16
Lab Sample ID:	JC19644-3	Date Received:	05/05/16
Matrix:	AQ - Ground Water	Percent Solids:	n/a
Method:	SW846 8270D BY SIM SW846 3510C		
Project:	BMSMC, Building 5 Area, PR		

	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	3M61130.D	1	05/06/16	LK	05/05/16	OP93662A	E3M2874
Run #2							

	Initial Volume	Final Volume
Run #1	900 ml	1.0 ml
Run #2		

CAS No.	Compound	Result	RL	MDL	Units	Q
91-20-3	Naphthalene	ND	0.11	0.033	ug/l	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
4165-60-0	Nitrobenzene-d5	66%		24-125%
321-60-8	2-Fluorobiphenyl	66%		19-127%
1718-51-0	Terphenyl-d14	78%		10-119%



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID:	RA-11 GWS	Date Sampled:	05/03/16
Lab Sample ID:	JC19644-3	Date Received:	05/05/16
Matrix:	AQ - Ground Water	Percent Solids:	n/a
Method:	SW846-8015C (DAI)		
Project:	BMSMC, Building 5 Area, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	GH104967.D	1	05/16/16	XPL	n/a	n/a	GGH5285
Run #2							

Low Molecular Alcohol List

CAS No.	Compound	Result	RL	MDL	Units	Q
64-17-5	Ethanol	ND	100	55	ug/l	
78-83-1	Isobutyl Alcohol	ND	100	36	ug/l	
67-63-0	Isopropyl Alcohol	ND	100	68	ug/l	
71-23-8	n-Propyl Alcohol	ND	100	43	ug/l	
71-36-3	n-Butyl Alcohol	ND	100	87	ug/l	
78-92-2	sec-Butyl Alcohol	ND	100	66	ug/l	
67-56-1	Methanol	ND	200	71	ug/l	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
111-27-3	Hexanol	105%		56-145%



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID:	RA-11 GWS	Date Sampled:	05/03/16
Lab Sample ID:	JC19644-3	Date Received:	05/05/16
Matrix:	AQ - Ground Water	Percent Solids:	n/a
Method:	SW846 8081B SW846 3510C		
Project:	BMSMC, Building 5 Area, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	4G67962.D	1	05/06/16	BP	05/05/16	OP93665	G4G1782
Run #2							

Run #	Initial Volume	Final Volume
Run #1	300 ml	2.0 ml
Run #2		

Pesticide TCL List

CAS No.	Compound	Result	RL	MDL	Units	Q
309-00-2	Aldrin	ND	0.0067	0.0040	ug/l	
319-84-6	alpha-BHC	ND	0.0067	0.0040	ug/l	
319-85-7	beta-BHC	ND	0.0067	0.0038	ug/l	
319-86-8	delta-BHC	ND	0.0067	0.0030	ug/l	
58-89-9	gamma-BHC (Lindane)	ND	0.0067	0.0019	ug/l	
5103-71-9	alpha-Chlordane	ND	0.0067	0.0031	ug/l	
5103-74-2	gamma-Chlordane	ND	0.0067	0.0031	ug/l	
60-57-1	Dieldrin	ND	0.0067	0.0024	ug/l	
72-54-8	4,4'-DDD	ND	0.0067	0.0025	ug/l	
72-55-9	4,4'-DDE	ND	0.0067	0.0041	ug/l	
50-29-3	4,4'-DDT	ND	0.0067	0.0033	ug/l	
72-20-8	Endrin	ND	0.0067	0.0034	ug/l	
1031-07-8	Endosulfan sulfate	ND	0.0067	0.0035	ug/l	
7421-93-4	Endrin aldehyde	ND	0.0067	0.0034	ug/l	
53494-70-5	Endrin ketone	ND	0.0067	0.0034	ug/l	
959-98-8	Endosulfan-I	ND	0.0067	0.0033	ug/l	
33213-65-9	Endosulfan-II	ND	0.0067	0.0029	ug/l	
76-44-8	Heptachlor	ND	0.0067	0.0025	ug/l	
1024-57-3	Heptachlor epoxide	ND	0.0067	0.0044	ug/l	
72-43-5	Methoxychlor	ND	0.013	0.0038	ug/l	
8001-35-2	Toxaphene	ND	0.17	0.12	ug/l	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
877-09-8	Tetrachloro-m-xylene	64%		26-132%
877-09-8	Tetrachloro-m-xylene	64%		26-132%
2051-24-3	Decachlorobiphenyl	97%		10-118%
2051-24-3	Decachlorobiphenyl	85%		10-118%



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 3

Client Sample ID:	RA-11 (10-11)	Date Sampled:	05/04/16
Lab Sample ID:	JC19644-4	Date Received:	05/05/16
Matrix:	SO - Soil	Percent Solids:	77.4
Method:	SW846 8270D SW846 3546		
Project:	BMSMC, Building 5 Area, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	3E83008.D	1	05/09/16	AN	05/06/16	OP93678	E3E3639
Run #2							

Run #	Initial Weight	Final Volume
Run #1	30.5 g	1.0 ml
Run #2		

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units	Q
95-57-8	2-Chlorophenol	ND	85	29	ug/kg	
59-50-7	4-Chloro-3-methyl phenol	ND	210	38	ug/kg	
120-83-2	2,4-Dichlorophenol	ND	210	42	ug/kg	
105-67-9	2,4-Dimethylphenol	ND	210	100	ug/kg	
51-28-5	2,4-Dinitrophenol	ND	210	170	ug/kg	
534-52-1	4,6-Dinitro-o-cresol	ND	210	62	ug/kg	
95-48-7	2-Methylphenol	ND	85	48	ug/kg	
	3&4-Methylphenol	ND	85	48	ug/kg	
88-75-5	2-Nitrophenol	ND	210	37	ug/kg	
100-02-7	4-Nitrophenol	ND	420	120	ug/kg	
87-86-5	Pentachlorophenol	ND	210	36	ug/kg	
108-95-2	Phenol	ND	85	33	ug/kg	
58-90-2	2,3,4,6-Tetrachlorophenol	ND	210	55	ug/kg	
95-95-4	2,4,5-Trichlorophenol	ND	210	61	ug/kg	
88-06-2	2,4,6-Trichlorophenol	ND	210	47	ug/kg	
83-32-9	Acenaphthene	ND	42	8.0	ug/kg	
208-96-8	Acenaphthylene	ND	42	5.8	ug/kg	
98-86-2	Acetophenone	ND	210	22	ug/kg	
120-12-7	Anthracene	ND	42	18	ug/kg	
1912-24-9	Atrazine	ND	85	14	ug/kg	
56-55-3	Benzo(a)anthracene	ND	42	6.5	ug/kg	
50-32-8	Benzo(a)pyrene	ND	42	10	ug/kg	
205-99-2	Benzo(b)fluoranthene	ND	42	8.9	ug/kg	
191-24-2	Benzo(g,h,i)perylene	ND	42	11	ug/kg	
207-08-9	Benzo(k)fluoranthene	ND	42	12	ug/kg	
101-55-3	4-Bromophenyl phenyl ether	ND	85	19	ug/kg	
85-68-7	Butyl benzyl phthalate	ND	85	24	ug/kg	
92-52-4	1,1'-Biphenyl	ND	85	12	ug/kg	
100-52-7	Benzaldehyde	ND	210	10	ug/kg	
91-58-7	2-Chloronaphthalene	ND	85	8.0	ug/kg	
106-47-8	4-Chloroaniline	ND	210	12	ug/kg	
86-74-8	Carbazole	ND	85	8.3	ug/kg	



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

Report of Analysis

Client Sample ID: RA-11 (10-11)
 Lab Sample ID: JC19644-4
 Matrix: SO - Soil
 Method: SW846 8270D SW846 3546
 Project: BSMC, Building 5 Area, PR

Date Sampled: 05/04/16
 Date Received: 05/05/16
 Percent Solids: 77.4

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units	Q
105-60-2	Caprolactam	ND	85	31	ug/kg	
218-01-9	Chrysene	ND	42	11	ug/kg	
111-91-1	bis(2-Chloroethoxy)methane	ND	85	8.9	ug/kg	
111-44-4	bis(2-Chloroethyl)ether	ND	85	24	ug/kg	
108-60-1	bis(2-Chloroisopropyl)ether	ND	85	17	ug/kg	
7005-72-3	4-Chlorophenyl phenyl ether	ND	85	11	ug/kg	
121-14-2	2,4-Dinitrotoluene	ND	42	15	ug/kg	
606-20-2	2,6-Dinitrotoluene	ND	42	13	ug/kg	
91-94-1	3,3'-Dichlorobenzidine	ND	85	57	ug/kg	
53-70-3	Dibenzo(a,h)anthracene	ND	42	16	ug/kg	
132-64-9	Dibenzofuran	ND	85	7.7	ug/kg	
84-74-2	Di-n-butyl phthalate	ND	85	26	ug/kg	
117-84-0	Di-n-octyl phthalate	ND	85	27	ug/kg	
84-66-2	Diethyl phthalate	ND	85	9.5	ug/kg	
131-11-3	Dimethyl phthalate	ND	85	8.7	ug/kg	
117-81-7	bis(2-Ethylhexyl)phthalate	ND	85	9.2	ug/kg	
206-44-0	Fluoranthene	ND	42	20	ug/kg	
86-73-7	Fluorene	ND	42	17	ug/kg	
118-74-1	Hexachlorobenzene	ND	85	11	ug/kg	
87-68-3	Hexachlorobutadiene	ND	42	17	ug/kg	
77-47-4	Hexachlorocyclopentadiene	ND	420	25	ug/kg	
67-72-1	Hexachloroethane	ND	210	32	ug/kg	
193-39-5	Indeno(1,2,3-cd)pyrene	ND	42	14	ug/kg	
78-59-1	Isophorone	ND	85	9.4	ug/kg	
90-12-0	1-Methylnaphthalene	ND	85	7.5	ug/kg	
91-57-6	2-Methylnaphthalene	ND	85	34	ug/kg	
88-74-4	2-Nitroaniline	ND	210	31	ug/kg	
99-09-2	3-Nitroaniline	ND	210	15	ug/kg	
100-01-6	4-Nitroaniline	ND	210	17	ug/kg	
98-95-3	Nitrobenzene	ND	85	19	ug/kg	
621-64-7	N-Nitroso-di-n-propylamine	ND	85	19	ug/kg	
86-30-6	N-Nitrosodiphenylamine	ND	210	24	ug/kg	
85-01-8	Phenanthrene	ND	42	10	ug/kg	
129-00-0	Pyrene	ND	42	7.4	ug/kg	
95-94-3	1,2,4,5-Tetrachlorobenzene	ND	210	9.5	ug/kg	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
367-12-4	2-Fluorophenol	70%		30-106%
4165-62-2	Phenol-d5	69%		30-106%

ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound



Report of Analysis

Client Sample ID: RA-11 (10-11)	
Lab Sample ID: JC19644-4	Date Sampled: 05/04/16
Matrix: SO - Soil	Date Received: 05/05/16
Method: SW846 8270D SW846 3546	Percent Solids: 77.4
Project: BSMC, Building 5 Area, PR	

ABN TCL Special List

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
118-79-6	2,4,6-Tribromophenol	83%		24-140%
4165-60-0	Nitrobenzene-d5	81%		26-122%
321-60-8	2-Fluorobiphenyl	89%		36-112%
1718-51-0	Terphenyl-d14	91%		36-132%



ND = Not detected MDL = Method Detection Limit
RL = Reporting Limit
E = Indicates value exceeds calibration range

J = Indicates an estimated value
B = Indicates analyte found in associated method blank
N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID:	RA-11 (10-11)	Date Sampled:	05/04/16
Lab Sample ID:	JC19644-4	Date Received:	05/05/16
Matrix:	SO - Soil	Percent Solids:	77.4
Method:	SW846 8270D BY SIM SW846 3546		
Project:	BMSMC, Building 5 Area, PR		

	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	4M65183.D	1	05/09/16	LK	05/06/16	OP93678A	E4M2910
Run #2							

	Initial Weight	Final Volume
Run #1	30.5 g	1.0 ml
Run #2		

CAS No.	Compound	Result	RL	MDL	Units	Q
123-91-1	1,4-Dioxane ^a	8.07	4.2	0.85	ug/kg	
91-20-3	Naphthalene	ND	4.2	0.52	ug/kg	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
4165-60-0	Nitrobenzene-d5	103%		15-138%
321-60-8	2-Fluorobiphenyl	74%		12-148%
1718-51-0	Terphenyl-d14	91%		10-157%

(a) Not accredited for this compound at the time of analysis, but all method requirements were followed.



ND = Not detected MDL = Method Detection Limit

RL = Reporting Limit

E = Indicates value exceeds calibration range

J = Indicates an estimated value

B = Indicates analyte found in associated method blank

N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID:	RA-11 (10-11)	Date Sampled:	05/04/16
Lab Sample ID:	JC19644-4	Date Received:	05/05/16
Matrix:	SO - Soil	Percent Solids:	77.4
Method:	SW846-8015C (DAD)		
Project:	BMSMC, Building 5 Area, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	GH105052.D	1	05/18/16	XPL	n/a	n/a	GGH5290
Run #2							

Run #	Initial Weight
Run #1	5.0 g
Run #2	

Low Molecular Alcohol List

CAS No.	Compound	Result	RL	MDL	Units	Q
64-17-5	Ethanol	ND	130	89	ug/kg	
78-83-1	Isobutyl Alcohol	ND	130	76	ug/kg	
67-63-0	Isopropyl Alcohol	ND	130	74	ug/kg	
71-23-8	n-Propyl Alcohol	ND	130	52	ug/kg	
71-36-3	n-Butyl Alcohol	ND	130	70	ug/kg	
78-92-2	sec-Butyl Alcohol	ND	130	69	ug/kg	
67-56-1	Methanol	ND	260	62	ug/kg	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
111-27-3	Hexanol	73%		52-141%
111-27-3	Hexanol	78%		52-141%



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID:	RA-11 (10-11)	Date Sampled:	05/04/16
Lab Sample ID:	JC19644-4	Date Received:	05/05/16
Matrix:	SO - Soil	Percent Solids:	77.4
Method:	SW846 8081B SW846 3546		
Project:	BMSMC, Building 5 Area, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	6G34906.D	1	05/10/16	RK	05/09/16	OP93736	G6G1007
Run #2							

Run #	Initial Weight	Final Volume
Run #1	16.3 g	10.0 ml
Run #2		

Pesticide TCL List

CAS No.	Compound	Result	RL	MDL	Units	Q
309-00-2	Aldrin	ND	0.79	0.71	ug/kg	
319-84-6	alpha-BHC	ND	0.79	0.53	ug/kg	
319-85-7	beta-BHC	ND	0.79	0.49	ug/kg	
319-86-8	delta-BHC	ND	0.79	0.31	ug/kg	
58-89-9	gamma-BHC (Lindane)	ND	0.79	0.36	ug/kg	
5103-71-9	alpha-Chlordane	ND	0.79	0.42	ug/kg	
5103-74-2	gamma-Chlordane	ND	0.79	0.60	ug/kg	
60-57-1	Dieldrin	ND	0.79	0.62	ug/kg	
72-54-8	4,4'-DDD	ND	0.79	0.29	ug/kg	
72-55-9	4,4'-DDE	ND	0.79	0.26	ug/kg	
50-29-3	4,4'-DDT	ND	0.79	0.30	ug/kg	
72-20-8	Endrin	ND	0.79	0.28	ug/kg	
1031-07-8	Endosulfan sulfate	ND	0.79	0.45	ug/kg	
7421-93-4	Endrin aldehyde	ND	0.79	0.59	ug/kg	
959-98-8	Endosulfan-I	ND	0.79	0.26	ug/kg	
33213-65-9	Endosulfan-II	ND	0.79	0.75	ug/kg	
76-44-8	Heptachlor	ND	0.79	0.65	ug/kg	
1024-57-3	Heptachlor epoxide	ND	0.79	0.33	ug/kg	
72-43-5	Methoxychlor	ND	1.6	0.44	ug/kg	
53494-70-5	Endrin ketone	ND	0.79	0.42	ug/kg	
8001-35-2	Toxaphene	ND	20	14	ug/kg	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
877-09-8	Tetrachloro-m-xylene	109%		24-136%
877-09-8	Tetrachloro-m-xylene	107%		24-136%
2051-24-3	Decachlorobiphenyl	106%		10-153%
2051-24-3	Decachlorobiphenyl	120%		10-153%



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 3

Client Sample ID:	RA-10 (5.5-6.5)	Date Sampled:	05/04/16
Lab Sample ID:	JC19644-5	Date Received:	05/05/16
Matrix:	SO - Soil	Percent Solids:	79.8
Method:	SW846 8270D SW846 3546		
Project:	BMSMC, Building 5 Area, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	3E83007.D	1	05/09/16	AN	05/06/16	OP93678	E3E3639
Run #2	3E83020.D	1	05/09/16	AN	05/06/16	OP93678	E3E3639

Run #	Initial Weight	Final Volume
Run #1	31.0 g	1.0 ml
Run #2	31.0 g	1.0 ml

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units	Q
95-57-8	2-Chlorophenol	ND	81	28	ug/kg	
59-50-7	4-Chloro-3-methyl phenol	ND	200	36	ug/kg	
120-83-2	2,4-Dichlorophenol	ND	200	40	ug/kg	
105-67-9	2,4-Dimethylphenol	ND	200	97	ug/kg	
51-28-5	2,4-Dinitrophenol	ND	200	160	ug/kg	
534-52-1	4,6-Dinitro-o-cresol	ND	200	59	ug/kg	
95-48-7	2-Methylphenol	ND	81	46	ug/kg	
	3&4-Methylphenol	ND	81	46	ug/kg	
88-75-5	2-Nitrophenol	ND	200	35	ug/kg	
100-02-7	4-Nitrophenol	ND	400	110	ug/kg	
87-86-5	Pentachlorophenol	ND	200	34	ug/kg	
108-95-2	Phenol	ND	81	31	ug/kg	
58-90-2	2,3,4,6-Tetrachlorophenol	ND	200	53	ug/kg	
95-95-4	2,4,5-Trichlorophenol	ND	200	58	ug/kg	
88-06-2	2,4,6-Trichlorophenol	ND	200	44	ug/kg	
83-32-9	Acenaphthene	ND	40	7.6	ug/kg	
208-96-8	Acenaphthylene	ND	40	5.6	ug/kg	
98-86-2	Acetophenone	ND	200	21	ug/kg	
120-12-7	Anthracene	ND	40	18	ug/kg	
1912-24-9	Atrazine	ND	81	13	ug/kg	
56-55-3	Benzo(a)anthracene	ND	40	6.2	ug/kg	
50-32-8	Benzo(a)pyrene	ND	40	9.9	ug/kg	
205-99-2	Benzo(b)fluoranthene	ND	40	8.5	ug/kg	
191-24-2	Benzo(g,h,i)perylene	ND	40	11	ug/kg	
207-08-9	Benzo(k)fluoranthene	ND	40	11	ug/kg	
101-55-3	4-Bromophenyl phenyl ether	ND	81	18	ug/kg	
85-68-7	Butyl benzyl phthalate	ND	81	23	ug/kg	
92-52-4	1,1'-Biphenyl	ND	81	11	ug/kg	
100-52-7	Benzaldehyde	ND	200	9.8	ug/kg	
91-58-7	2-Chloronaphthalene	ND	81	7.6	ug/kg	
106-47-8	4-Chloroaniline	ND	200	12	ug/kg	
86-74-8	Carbazole	ND	81	7.9	ug/kg	



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

Report of Analysis

Client Sample ID: RA-10 (5.5-6.5)
 Lab Sample ID: JC19644-5
 Matrix: SO - Soil
 Method: SW846 8270D SW846 3546
 Project: BSMC, Building 5 Area, PR

Date Sampled: 05/04/16
 Date Received: 05/05/16
 Percent Solids: 79.8

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units	Q
105-60-2	Caprolactam	ND	81	29	ug/kg	
218-01-9	Chrysene	ND	40	10	ug/kg	
111-91-1	bis(2-Chloroethoxy)methane	ND	81	8.5	ug/kg	
111-44-4	bis(2-Chloroethyl)ether	ND	81	23	ug/kg	
108-60-1	bis(2-Chloroisopropyl)ether	ND	81	16	ug/kg	
7005-72-3	4-Chlorophenyl phenyl ether	ND	81	11	ug/kg	
121-14-2	2,4-Dinitrotoluene	ND	40	15	ug/kg	
606-20-2	2,6-Dinitrotoluene	ND	40	12	ug/kg	
91-94-1	3,3'-Dichlorobenzidine	ND	81	55	ug/kg	
53-70-3	Dibenzo(a,h)anthracene	ND	40	15	ug/kg	
132-64-9	Dibenzofuran	ND	81	7.3	ug/kg	
84-74-2	Di-n-butyl phthalate	ND	81	25	ug/kg	
117-84-0	Di-n-octyl phthalate	ND	81	26	ug/kg	
84-66-2	Diethyl phthalate	ND	81	9.1	ug/kg	
131-11-3	Dimethyl phthalate	ND	81	8.3	ug/kg	
117-81-7	bis(2-Ethylhexyl)phthalate	ND	81	8.8	ug/kg	
206-44-0	Fluoranthene	ND	40	19	ug/kg	
86-73-7	Fluorene	ND	40	16	ug/kg	
118-74-1	Hexachlorobenzene	ND	81	10	ug/kg	
87-68-3	Hexachlorobutadiene	ND	40	16	ug/kg	
77-47-4	Hexachlorocyclopentadiene	ND	400	24	ug/kg	
67-72-1	Hexachloroethane	ND	200	31	ug/kg	
193-39-5	Indeno(1,2,3-cd)pyrene	ND	40	14	ug/kg	
78-59-1	Isophorone	ND	81	9.0	ug/kg	
90-12-0	1-Methylnaphthalene	ND	81	7.2	ug/kg	
91-57-6	2-Methylnaphthalene	ND	81	33	ug/kg	
88-74-4	2-Nitroaniline	ND	200	29	ug/kg	
99-09-2	3-Nitroaniline	ND	200	14	ug/kg	
100-01-6	4-Nitroaniline	ND	200	16	ug/kg	
98-95-3	Nitrobenzene	ND	81	18	ug/kg	
621-64-7	N-Nitroso-di-n-propylamine	ND	81	18	ug/kg	
86-30-6	N-Nitrosodiphenylamine	ND	200	23	ug/kg	
85-01-8	Phenanthrene	ND	40	9.7	ug/kg	
129-00-0	Pyrene	ND	40	7.1	ug/kg	
95-94-3	1,2,4,5-Tetrachlorobenzene	ND	200	9.1	ug/kg	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
367-12-4	2-Fluorophenol	70%	66%	30-106%
4165-62-2	Phenol-d5	68%	66%	30-106%

ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound



Report of Analysis

Client Sample ID:	RA-10 (5.5-6.5)	Date Sampled:	05/04/16
Lab Sample ID:	JC19644-5	Date Received:	05/05/16
Matrix:	SO - Soil	Percent Solids:	79.8
Method:	SW846 8270D SW846 3546		
Project:	BMSMC, Building 5 Area, PR		

ABN TCL Special List

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
118-79-6	2,4,6-Tribromophenol	64%	62%	24-140%
4165-60-0	Nitrobenzene-d5	80%	79%	26-122%
321-60-8	2-Fluorobiphenyl	80%	79%	36-112%
1718-51-0	Terphenyl-d14	84%	90%	36-132%

(a) Confirmation run for internal standard areas.



ND = Not detected MDL = Method Detection Limit
RL = Reporting Limit
E = Indicates value exceeds calibration range

J = Indicates an estimated value
B = Indicates analyte found in associated method blank
N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID: RA-10 (5.5-6.5)
 Lab Sample ID: JC19644-5
 Matrix: SO - Soil
 Method: SW846 8270D BY SIM SW846 3546
 Project: BSMC, Building 5 Area, PR

Date Sampled: 05/04/16
 Date Received: 05/05/16
 Percent Solids: 79.8

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	4M65184.D	1	05/09/16	LK	05/06/16	OP93678A	E4M2910
Run #2							

Run #	Initial Weight	Final Volume
Run #1	31.0 g	1.0 ml
Run #2		

CAS No.	Compound	Result	RL	MDL	Units	Q
123-91-1	1,4-Dioxane ^a	136	4.0	0.81	ug/kg	
91-20-3	Naphthalene	ND	4.0	0.49	ug/kg	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
4165-60-0	Nitrobenzene-d5	91%		15-138%
321-60-8	2-Fluorobiphenyl	79%		12-148%
1718-51-0	Terphenyl-d14	90%		10-157%

(a) Not accredited for this compound at the time of analysis, but all method requirements were followed.



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID:	RA-10 (5.5-6.5)	Date Sampled:	05/04/16
Lab Sample ID:	JC19644-5	Date Received:	05/05/16
Matrix:	SO - Soil	Percent Solids:	79.8
Method:	SW846-8015C (DAI)		
Project:	BMSMC, Building 5 Area, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	GH105053.D	1	05/18/16	XPL	n/a	n/a	GGH5290
Run #2							

Run #	Initial Weight
Run #1	5.0 g
Run #2	

Low Molecular Alcohol List

CAS No.	Compound	Result	RL	MDL	Units	Q
64-17-5	Ethanol	ND	130	86	ug/kg	
78-83-1	Isobutyl Alcohol	ND	130	74	ug/kg	
67-63-0	Isopropyl Alcohol	ND	130	72	ug/kg	
71-23-8	n-Propyl Alcohol	ND	130	50	ug/kg	
71-36-3	n-Butyl Alcohol	ND	130	68	ug/kg	
78-92-2	sec-Butyl Alcohol	ND	130	67	ug/kg	
67-56-1	Methanol	ND	250	60	ug/kg	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
111-27-3	Hexanol	84%		52-141%
111-27-3	Hexanol	88%		52-141%



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID: RA-10 (5.5-6.5)
 Lab Sample ID: JC19644-5
 Matrix: SO - Soil
 Method: SW846 8081B SW846 3546
 Project: BSMC, Building 5 Area, PR

Date Sampled: 05/04/16
 Date Received: 05/05/16
 Percent Solids: 79.8

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	6G34907.D	1	05/10/16	RK	05/09/16	OP93736	G6G1007
Run #2							

Run #	Initial Weight	Final Volume
Run #1	15.1 g	10.0 ml
Run #2		

Pesticide TCL List

CAS No.	Compound	Result	RL	MDL	Units	Q
309-00-2	Aldrin	ND	0.83	0.74	ug/kg	
319-84-6	alpha-BHC	ND	0.83	0.55	ug/kg	
319-85-7	beta-BHC	ND	0.83	0.51	ug/kg	
319-86-8	delta-BHC	ND	0.83	0.33	ug/kg	
58-89-9	gamma-BHC (Lindane)	ND	0.83	0.38	ug/kg	
5103-71-9	alpha-Chlordane	ND	0.83	0.44	ug/kg	
5103-74-2	gamma-Chlordane	ND	0.83	0.63	ug/kg	
60-57-1	Dieldrin	ND	0.83	0.65	ug/kg	
72-54-8	4,4'-DDD	ND	0.83	0.31	ug/kg	
72-55-9	4,4'-DDE	ND	0.83	0.28	ug/kg	
50-29-3	4,4'-DDT	ND	0.83	0.32	ug/kg	
72-20-8	Endrin	ND	0.83	0.29	ug/kg	
1031-07-8	Endosulfan sulfate	ND	0.83	0.47	ug/kg	
7421-93-4	Endrin aldehyde	ND	0.83	0.62	ug/kg	
959-98-8	Endosulfan-I	ND	0.83	0.27	ug/kg	
33213-65-9	Endosulfan-II	ND	0.83	0.79	ug/kg	
76-44-8	Heptachlor	ND	0.83	0.68	ug/kg	
1024-57-3	Heptachlor epoxide	ND	0.83	0.34	ug/kg	
72-43-5	Methoxychlor	ND	1.7	0.46	ug/kg	
53494-70-5	Endrin ketone	ND	0.83	0.44	ug/kg	
8001-35-2	Toxaphene	ND	21	14	ug/kg	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
877-09-8	Tetrachloro-m-xylene	104%		24-136%
877-09-8	Tetrachloro-m-xylene	107%		24-136%
2051-24-3	Decachlorobiphenyl	109%		10-153%
2051-24-3	Decachlorobiphenyl	119%		10-153%



ND = Not detected MDL = Method Detection Limit
 RL = Reporting Limit
 E = Indicates value exceeds calibration range

J = Indicates an estimated value
 B = Indicates analyte found in associated method blank
 N = Indicates presumptive evidence of a compound

SGS

GW
SO- Foreign
ACCUTEST - NJ

CHAIN OF CUSTODY

SGS Accutest - Dayton
2235 Route 130, Dayton, NJ 08510
TEL: 732-329-4200 FAX: 732-329-3499/3480
www.accutest.com

PAGE 1 OF 1

510

801219536007

SGS Account Card #

SGS Account Card # JC19644

Company Name Anderson Mulholland Assoc. Inc.		Project Name BMS Release Assessment		Matrix Codes	
Address 2700 Westchester		City Humacao State PR		DW - Drinking Water GW - Ground Water WW - Water SW - Surface Water SD - Soil SL - Sludge SED - Sediment CI - Oil LIQ - Other Liquid AIR - Air SOL - Other Solid WP - Waste PB Field Blank EB - Equipment Blank RB - Reuse Blank TB - Trip Blank	
City Purchase NY		Billing Information (if different from Report to)			
Phone Contact Terry Taylor		Company Name			
Phone # 914-251-0400		Street Address			
Sample(s) Name(s) N. Rivera, T. Taylor, B. Lindstrom		City			
Phone #		State			
Project Manager		Zip			
Client Purchase Order #		Admission			
Collection		Number of preserved bottles			
Field ID / Point of Collection	Date	Time	Sampled by	Matrix	# of bottles
1 RA-16 GWD	5/3/16	1100	TT GW	7	3
2 S-425	5/3/16	1535	NR GW	7	3
3 RA-11 GWS	5/3/16	1645	TT GW	7	3
4 RA-11 (10-11)	5/4/16	1100	TT SO	2	2
5 RA-10 (5.5-6.5)	5/4/16	1235	NR SO	2	2
INITIAL ASSESSMENT 2/33 LABEL VERIFICATION 710					
Turnaround Time / Business days: _____ Data Deliverable Information: _____ Comments / Special Instructions: Add to Report: 1-methyl naphthalene with SVOC by Method 8270D					
Emergency & Rush T/A rate applies via Letter _____ Sample Custody must be maintained before each time sample change processing, including courier delivery.					
1	5/3/16 1500	1	FEDEX EX	2	FEDEX
2		2		3	
3		3		4	
4		4		5	
5		5			

JC19644: Chain of Custody

Page 1 of 2

SGS

40 of 1305
ACCUTEST
JC19644

EXECUTIVE NARRATIVE

SDG No: **JC19644** Laboratory: **Accutest, New Jersey**
Analysis: **SW846-8270D** Number of Samples: **5**
Location: **BMSMC, Building 5 Area**
Humacao, PR

SUMMARY: Five (5) samples were analyzed for the ABN TCL list following method SW846-8270D; Naphthalene and 1,4-Dioxane were also analyzed by SW846-8270D using the selective ion monitoring (SIM) technique. The sample results were assessed according to USEPA data validation guidance documents in the following order of precedence: EPA Hazardous Waste Support Section, SOP HW-35A, July 2015 –Revision 0. *Semivolatile Data Validation*. The QC criteria and data validation actions listed on the data review worksheets are from the primary guidance document, unless otherwise noted.

Results are valid and can be used for decision making purposes.

Critical issues: **None**
Major: **None**
Minor: **None**

Critical findings: **None**
Major findings: **None**
Minor findings:

1 Initial and continuing calibration verifications meet the method and guidance document required performance criteria except the cases described above. Results qualified as estimated (J), (UJ) for non-detects.

No closing calibration verification included in data package. No action taken, professional judgment.

No qualification performed on analytes with % difference in the continue calibration verification outside the method performance criteria but within the validation guidelines criteria, + 40 %. No action taken.

QC samples for this data package were analyzed in the instrument GCM6P (Scan). Analytes not meeting continue calibration verification criteria were not qualified.

2. Surrogate (acid fraction) not recovered in sample JC19644-2 due to dilution. No action taken.


3. Hexachlorocyclopentadiene MS/MSD results were outside the laboratory control limit. No action was taken, MS/MSD results apply to unspiked sample. Unspiked sample was from a different project.

4. Internal standards perylene-d12 and chrysene-d12 outside the laboratory control limit in sample JC19644-5. No action taken.

COMMENTS: Results are valid and can be used for decision making purposes.

Reviewers Name: Rafael Infante
Chemist License 1888

Signature:

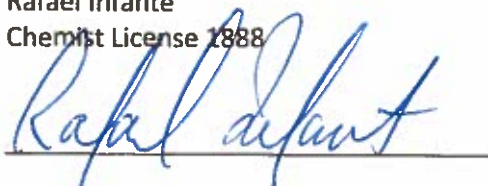

Date: May 30, 2016

SAMPLE ORGANIC DATA SAMPLE SUMMARY

Sample ID: JC19644-1

Sample location: BMSMC Building 5 Area

Sampling date: 5/3/2016

Matrix: Groundwater

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
2-Chlorophenol	5.6	ug/l	1	-	U	Yes
4-Chloro-3-methyl phenol	5.6	ug/l	1	-	U	Yes
2,4-Dichlorophenol	2.2	ug/l	1	-	U	Yes
2,4-Dimethylphenol	5.6	ug/l	1	-	U	Yes
2,4-Dinitrophenol	11	ug/l	1	-	U	Yes
4,6-Dinitro-o-cresol	5.6	ug/l	1	-	U	Yes
2-Methylphenol	2.2	ug/l	1	-	U	Yes
3&4-Methylphenol	2.2	ug/l	1	-	U	Yes
2-Nitrophenol	5.6	ug/l	1	-	U	Yes
4-Nitrophenol	11	ug/l	1	-	U	Yes
Pentachlorophenol	5.6	ug/l	1	-	U	Yes
Phenol	2.2	ug/l	1	-	U	Yes
2,3,4,6-Tetrachlorophenol	5.6	ug/l	1	-	U	Yes
2,4,5-Trichlorophenol	5.6	ug/l	1	-	U	Yes
2,4,6-Trichlorophenol	5.6	ug/l	1	-	U	Yes
Acenaphthene	1.0	ug/l	1	-	U	Yes
Acenaphthylene	1.0	ug/l	1	-	U	Yes
Acetophenone	2.2	ug/l	1	-	U	Yes
Anthracene	1.0	ug/l	1	-	U	Yes
Atrazine	2.2	ug/l	1	-	U	Yes
Benzaldehyde	5.6	ug/l	1	-	U	Yes
Benzo(a)anthracene	1.0	ug/l	1	-	U	Yes
Benzo(a)pyrene	1.0	ug/l	1	-	U	Yes
Benzo(b)fluoranthene	1.0	ug/l	1	-	U	Yes
Benzo(g,h,i)perylene	1.0	ug/l	1	-	U	Yes
Benzo(k)fluoranthene	1.0	ug/l	1	-	U	Yes
4-Bromophenyl phenyl ether	2.2	ug/l	1	-	U	Yes
Butyl benzyl phthalate	2.2	ug/l	1	-	U	Yes
1,1'-Biphenyl	1.0	ug/l	1	-	U	Yes
2-Chloronaphthalene	2.2	ug/l	1	-	U	Yes
4-Chloroaniline	5.6	ug/l	1	-	U	Yes
Carbazole	1.0	ug/l	1	-	U	Yes
Caprolactam	2.2	ug/l	1	-	U	Yes
Chrysene	1.0	ug/l	1	-	U	Yes
bis(2-Chloroethoxy)methane	2.2	ug/l	1	-	U	Yes
bis(2-Chloroethyl)ether	2.2	ug/l	1	-	U	Yes

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
bis(2-Chloroisopropyl)ether	2.2	ug/l	1	-	U	Yes
4-Chlorophenyl phenyl ether	2.2	ug/l	1	-	U	Yes
2,4-Dinitrotoluene	1.0	ug/l	1	-	U	Yes
2,6-Dinitrotoluene	1.0	ug/l	1	-	U	Yes
3,3'-Dichlorobenzidine	2.2	ug/l	1	-	U	Yes
1,4-Dioxane	26.7	ug/l	1	-	-	Yes
Dibenzo(a,h)anthracene	1.0	ug/l	1	-	U	Yes
Dibenzofuran	5.6	ug/l	1	-	U	Yes
Di-n-butyl phthalate	2.2	ug/l	1	-	U	Yes
Di-n-octyl phthalate	2.2	ug/l	1	-	U	Yes
Diethyl phthalate	2.2	ug/l	1	-	U	Yes
Dimethyl phthalate	2.2	ug/l	1	-	U	Yes
bis(2-Ethylhexyl)phthalate	2.2	ug/l	1	-	U	Yes
Fluoranthene	1.0	ug/l	1	-	U	Yes
Fluorene	1.0	ug/l	1	-	U	Yes
Hexachlorobenzene	1.0	ug/l	1	-	U	Yes
Hexachlorobutadiene	1.0	ug/l	1	-	U	Yes
Hexachlorocyclopentadiene	11	ug/l	1	-	U	Yes
Hexachloroethane	2.2	ug/l	1	-	U	Yes
Indeno(1,2,3-cd)pyrene	1.0	ug/l	1	-	U	Yes
Isophorone	2.2	ug/l	1	-	U	Yes
1-Methylnaphthalene	1.0	ug/l	1	-	U	Yes
2-Methylnaphthalene	1.0	ug/l	1	-	U	Yes
2-Nitroaniline	5.6	ug/l	1	-	UJ	Yes
3-Nitroaniline	5.6	ug/l	1	-	U	Yes
4-Nitroaniline	5.6	ug/l	1	-	U	Yes
Nitrobenzene	2.2	ug/l	1	-	U	Yes
N-Nitroso-di-n-propylamine	2.2	ug/l	1	-	UJ	Yes
Nitrosodiphenylamine	5.6	ug/l	1	-	U	Yes
Phenanthrene	1.0	ug/l	1	-	U	Yes
Pyrene	1.0	ug/l	1	-	U	Yes
1,2,4,5-Tetrachlorobenzene	2.2	ug/l	1	-	U	Yes

METHOD: 8270D (SIM)

Naphthalene	0.11	ug/l	1	-	U	Yes
-------------	------	------	---	---	---	-----

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
--------------	--------	-------	-----------------	----------	------------	------------

Sample ID: JC19644-2

Sample location: BMSMC Building 5 Area

Sampling date: 5/3/2016

Matrix: Groundwater

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
2-Chlorophenol	5.0	ug/l	1	-	U	Yes
4-Chloro-3-methyl phenol	5.0	ug/l	1	-	U	Yes
2,4-Dichlorophenol	2.0	ug/l	1	-	U	Yes
2,4-Dimethylphenol	5.0	ug/l	1	-	U	Yes
2,4-Dinitrophenol	10	ug/l	1	-	U	Yes
4,6-Dinitro-o-cresol	5.0	ug/l	1	-	U	Yes
2-Methylphenol	2.0	ug/l	1	-	U	Yes
3&4-Methylphenol	2.0	ug/l	1	-	U	Yes
2-Nitrophenol	5.0	ug/l	1	-	U	Yes
4-Nitrophenol	10	ug/l	1	-	U	Yes
Pentachlorophenol	5.0	ug/l	1	-	U	Yes
Phenol	2.0	ug/l	1	-	U	Yes
2,3,4,6-Tetrachlorophenol	5.0	ug/l	1	-	U	Yes
2,4,5-Trichlorophenol	5.0	ug/l	1	-	U	Yes
2,4,6-Trichlorophenol	5.0	ug/l	1	-	U	Yes
Acenaphthene	1.0	ug/l	1	-	U	Yes
Acenaphthylene	1.0	ug/l	1	-	U	Yes
Acetophenone	2.0	ug/l	1	-	U	Yes
Anthracene	1.0	ug/l	1	-	U	Yes
Atrazine	2.0	ug/l	1	-	U	Yes
Benzaldehyde	0.84	ug/l	1	J	UJ	Yes
Benzo(a)anthracene	1.0	ug/l	1	-	U	Yes
Benzo(a)pyrene	1.0	ug/l	1	-	U	Yes
Benzo(b)fluoranthene	1.0	ug/l	1	-	U	Yes
Benzo(g,h,i)perylene	1.0	ug/l	1	-	U	Yes
Benzo(k)fluoranthene	1.0	ug/l	1	-	U	Yes
4-Bromophenyl phenyl ether	2.0	ug/l	1	-	U	Yes
Butyl benzyl phthalate	2.0	ug/l	1	-	U	Yes
1,1'-Biphenyl	1.0	ug/l	1	-	U	Yes
2-Chloronaphthalene	2.0	ug/l	1	-	U	Yes
4-Chloroaniline	5.0	ug/l	1	-	U	Yes
Carbazole	1.0	ug/l	1	-	U	Yes
Caprolactam	2.0	ug/l	1	-	U	Yes
Chrysene	1.0	ug/l	1	-	U	Yes
bis(2-Chloroethoxy)methane	2.0	ug/l	1	-	U	Yes

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
bis(2-Chloroethyl)ether	2.0	ug/l	1	-	U	Yes
bis(2-Chloroisopropyl)ether	2.0	ug/l	1	-	U	Yes
4-Chlorophenyl phenyl ether	2.0	ug/l	1	-	U	Yes
2,4-Dinitrotoluene	1.0	ug/l	1	-	U	Yes
2,6-Dinitrotoluene	1.0	ug/l	1	-	U	Yes
3,3'-Dichlorobenzidine	2.0	ug/l	1	-	U	Yes
1,4-Dioxane	3890	ug/l	50	-	-	Yes
Dibenzo(a,h)anthracene	1.0	ug/l	1	-	U	Yes
Dibenzofuran	5.0	ug/l	1	-	U	Yes
Di-n-butyl phthalate	2.0	ug/l	1	-	U	Yes
Di-n-octyl phthalate	2.0	ug/l	1	-	U	Yes
Diethyl phthalate	2.0	ug/l	1	-	U	Yes
Dimethyl phthalate	2.0	ug/l	1	-	U	Yes
bis(2-Ethylhexyl)phthalate	2.0	ug/l	1	-	U	Yes
Fluoranthene	1.0	ug/l	1	-	U	Yes
Fluorene	1.0	ug/l	1	-	U	Yes
Hexachlorobenzene	1.0	ug/l	1	-	U	Yes
Hexachlorobutadiene	1.0	ug/l	1	-	U	Yes
Hexachlorocyclopentadiene	11	ug/l	1	-	U	Yes
Hexachloroethane	2.0	ug/l	1	-	U	Yes
Indeno(1,2,3-cd)pyrene	1.0	ug/l	1	-	U	Yes
Isophorone	2.0	ug/l	1	-	U	Yes
1-Methylnaphthalene	1.0	ug/l	1	-	U	Yes
2-Methylnaphthalene	1.0	ug/l	1	-	U	Yes
2-Nitroaniline	5.0	ug/l	1	-	UJ	Yes
3-Nitroaniline	5.0	ug/l	1	-	U	Yes
4-Nitroaniline	5.0	ug/l	1	-	U	Yes
Nitrobenzene	2.0	ug/l	1	-	U	Yes
N-Nitroso-di-n-propylamine	2.0	ug/l	1	-	UJ	Yes
Nitrosodiphenylamine	5.0	ug/l	1	-	U	Yes
Phenanthrene	1.0	ug/l	1	-	U	Yes
Pyrene	1.0	ug/l	1	-	U	Yes
1,2,4,5-Tetrachlorobenzene	2.0	ug/l	1	-	U	Yes

METHOD: 8270D (SIM)

Naphthalene	0.10	ug/l	1	-	U	Yes
-------------	------	------	---	---	---	-----

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
--------------	--------	-------	-----------------	----------	------------	------------

Sample ID: JC19644-3

Sample location: BMSMC Building 5 Area

Sampling date: 5/3/2016

Matrix: Groundwater

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
2-Chlorophenol	5.6	ug/l	1	-	U	Yes
4-Chloro-3-methyl phenol	5.6	ug/l	1	-	U	Yes
2,4-Dichlorophenol	2.2	ug/l	1	-	U	Yes
2,4-Dimethylphenol	5.6	ug/l	1	-	U	Yes
2,4-Dinitrophenol	11	ug/l	1	-	U	Yes
4,6-Dinitro-o-cresol	5.6	ug/l	1	-	U	Yes
2-Methylphenol	2.2	ug/l	1	-	U	Yes
3&4-Methylphenol	2.2	ug/l	1	-	U	Yes
2-Nitrophenol	5.6	ug/l	1	-	U	Yes
4-Nitrophenol	11	ug/l	1	-	U	Yes
Pentachlorophenol	5.6	ug/l	1	-	U	Yes
Phenol	2.2	ug/l	1	-	U	Yes
2,3,4,6-Tetrachlorophenol	5.6	ug/l	1	-	U	Yes
2,4,5-Trichlorophenol	5.6	ug/l	1	-	U	Yes
2,4,6-Trichlorophenol	5.6	ug/l	1	-	U	Yes
Acenaphthene	1.0	ug/l	1	-	U	Yes
Acenaphthylene	1.0	ug/l	1	-	U	Yes
Acetophenone	2.2	ug/l	1	-	U	Yes
Anthracene	1.0	ug/l	1	-	U	Yes
Atrazine	2.2	ug/l	1	-	U	Yes
Benzaldehyde	5.6	ug/l	1	-	U	Yes
Benzo(a)anthracene	1.0	ug/l	1	-	U	Yes
Benzo(a)pyrene	1.0	ug/l	1	-	U	Yes
Benzo(b)fluoranthene	1.0	ug/l	1	-	U	Yes
Benzo(g,h,i)perylene	1.0	ug/l	1	-	U	Yes
Benzo(k)fluoranthene	1.0	ug/l	1	-	U	Yes
4-Bromophenyl phenyl ether	2.2	ug/l	1	-	U	Yes
Butyl benzyl phthalate	2.2	ug/l	1	-	U	Yes
1,1'-Biphenyl	1.0	ug/l	1	-	U	Yes
2-Chloronaphthalene	2.2	ug/l	1	-	U	Yes
4-Chloroaniline	5.6	ug/l	1	-	U	Yes
Carbazole	1.0	ug/l	1	-	U	Yes
Caprolactam	2.2	ug/l	1	-	U	Yes
Chrysene	1.0	ug/l	1	-	U	Yes
bis(2-Chloroethoxy)methane	2.2	ug/l	1	-	U	Yes
bis(2-Chloroethyl)ether	2.2	ug/l	1	-	U	Yes

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
bis(2-Chloroisopropyl)ether	2.2	ug/l	1	-	U	Yes
4-Chlorophenyl phenyl ether	2.2	ug/l	1	-	U	Yes
2,4-Dinitrotoluene	1.0	ug/l	1	-	U	Yes
2,6-Dinitrotoluene	1.0	ug/l	1	-	U	Yes
3,3'-Dichlorobenzidine	2.2	ug/l	1	-	U	Yes
1,4-Dioxane	67.1	ug/l	1	-	-	Yes
Dibenzo(a,h)anthracene	1.0	ug/l	1	-	U	Yes
Dibenzofuran	5.6	ug/l	1	-	U	Yes
Di-n-butyl phthalate	2.2	ug/l	1	-	U	Yes
Di-n-octyl phthalate	2.2	ug/l	1	-	U	Yes
Diethyl phthalate	2.2	ug/l	1	-	U	Yes
Dimethyl phthalate	2.2	ug/l	1	-	U	Yes
bis(2-Ethylhexyl)phthalate	2.2	ug/l	1	-	U	Yes
Fluoranthene	1.0	ug/l	1	-	U	Yes
Fluorene	1.0	ug/l	1	-	U	Yes
Hexachlorobenzene	1.0	ug/l	1	-	U	Yes
Hexachlorobutadiene	1.0	ug/l	1	-	U	Yes
Hexachlorocyclopentadiene	11	ug/l	1	-	U	Yes
Hexachloroethane	2.2	ug/l	1	-	U	Yes
Indeno(1,2,3-cd)pyrene	1.0	ug/l	1	-	U	Yes
Isophorone	2.2	ug/l	1	-	U	Yes
1-Methylnaphthalene	1.0	ug/l	1	-	U	Yes
2-Methylnaphthalene	1.0	ug/l	1	-	U	Yes
2-Nitroaniline	5.6	ug/l	1	-	UJ	Yes
3-Nitroaniline	5.6	ug/l	1	-	U	Yes
4-Nitroaniline	5.6	ug/l	1	-	U	Yes
Nitrobenzene	2.2	ug/l	1	-	U	Yes
N-Nitroso-di-n-propylamine	2.2	ug/l	1	-	UJ	Yes
Nitrosodiphenylamine	5.6	ug/l	1	-	U	Yes
Phenanthrene	1.0	ug/l	1	-	U	Yes
Pyrene	1.0	ug/l	1	-	U	Yes
1,2,4,5-Tetrachlorobenzene	2.2	ug/l	1	-	U	Yes

METHOD: 8270D (SIM)

Naphthalene	0.11	ug/l	1	-	U	Yes
-------------	------	------	---	---	---	-----

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
--------------	--------	-------	-----------------	----------	------------	------------

Sample ID: JC19644-4

Sample location: BMSMC Building 5 Area

Sampling date: 5/4/2016

Matrix: Soil

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
2-Chlorophenol	85	ug/kg	1	-	U	Yes
4-Chloro-3-methyl phenol	210	ug/kg	1	-	U	Yes
2,4-Dichlorophenol	210	ug/kg	1	-	U	Yes
2,4-Dimethylphenol	210	ug/kg	1	-	U	Yes
2,4-Dinitrophenol	210	ug/kg	1	-	U	Yes
4,6-Dinitro-o-cresol	210	ug/kg	1	-	U	Yes
2-Methylphenol	85	ug/kg	1	-	U	Yes
3&4-Methylphenol	85	ug/kg	1	-	U	Yes
2-Nitrophenol	210	ug/kg	1	-	U	Yes
4-Nitrophenol	420	ug/kg	1	-	U	Yes
Pentachlorophenol	210	ug/kg	1	-	U	Yes
Phenol	85	ug/kg	1	-	U	Yes
2,3,4,6-Tetrachlorophenol	210	ug/kg	1	-	U	Yes
2,4,5-Trichlorophenol	210	ug/kg	1	-	U	Yes
2,4,6-Trichlorophenol	210	ug/kg	1	-	U	Yes
Acenaphthene	42	ug/kg	1	-	U	Yes
Acenaphthylene	42	ug/kg	1	-	U	Yes
Acetophenone	210	ug/kg	1	-	U	Yes
Anthracene	42	ug/kg	1	-	U	Yes
Atrazine	85	ug/kg	1	-	U	Yes
Benzo(a)anthracene	42	ug/kg	1	-	U	Yes
Benzo(a)pyrene	42	ug/kg	1	-	U	Yes
Benzo(b)fluoranthene	42	ug/kg	1	-	U	Yes
Benzo(g,h,i)perylene	42	ug/kg	1	-	U	Yes
Benzo(k)fluoranthene	42	ug/kg	1	-	U	Yes
4-Bromophenyl phenyl ether	85	ug/kg	1	-	U	Yes
Butyl benzyl phthalate	85	ug/kg	1	-	U	Yes
1,1'-Biphenyl	85	ug/kg	1	-	U	Yes
Benzaldehyde	210	ug/kg	1	-	U	Yes
2-Chloronaphthalene	85	ug/kg	1	-	U	Yes
4-Chloroaniline	210	ug/kg	1	-	U	Yes
Carbazole	85	ug/kg	1	-	U	Yes
Caprolactam	85	ug/kg	1	-	U	Yes
Chrysene	42	ug/kg	1	-	U	Yes
bis(2-Chloroethoxy)methane	85	ug/kg	1	-	U	Yes
bis(2-Chloroethyl)ether	85	ug/kg	1	-	U	Yes

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
bis(2-Chloroisopropyl)ether	85	ug/kg	1	-	U	Yes
4-Chlorophenyl phenyl ether	85	ug/kg	1	-	U	Yes
2,4-Dinitrotoluene	42	ug/kg	1	-	U	Yes
2,6-Dinitrotoluene	42	ug/kg	1	-	U	Yes
3,3'-Dichlorobenzidine	85	ug/kg	1	-	U	Yes
Dibenzo(a,h)anthracene	42	ug/kg	1	-	U	Yes
Dibenzofuran	85	ug/kg	1	-	U	Yes
Di-n-butyl phthalate	85	ug/kg	1	-	U	Yes
Di-n-octyl phthalate	85	ug/kg	1	-	U	Yes
Diethyl phthalate	85	ug/kg	1	-	U	Yes
Dimethyl phthalate	85	ug/kg	1	-	U	Yes
bis(2-Ethylhexyl)phthalate	85	ug/kg	1	-	U	Yes
Fluoranthene	42	ug/kg	1	-	U	Yes
Fluorene	42	ug/kg	1	-	U	Yes
Hexachlorobenzene	85	ug/kg	1	-	U	Yes
Hexachlorobutadiene	42	ug/kg	1	-	UJ	Yes
Hexachlorocyclopentadiene	420	ug/kg	1	-	U	Yes
Hexachloroethane	210	ug/kg	1	-	U	Yes
Indeno(1,2,3-cd)pyrene	42	ug/kg	1	-	U	Yes
Isophorone	85	ug/kg	1	-	U	Yes
1-Methylnaphthalene	85	ug/kg	1	-	U	Yes
2-Methylnaphthalene	85	ug/kg	1	-	U	Yes
2-Nitroaniline	210	ug/kg	1	-	U	Yes
3-Nitroaniline	210	ug/kg	1	-	U	Yes
4-Nitroaniline	210	ug/kg	1	-	U	Yes
Nitrobenzene	85	ug/kg	1	-	U	Yes
N-Nitroso-di-n-propylamine	85	ug/kg	1	-	U	Yes
Nitrosodiphenylamine	210	ug/kg	1	-	U	Yes
Phenanthrene	42	ug/kg	1	-	U	Yes
Pyrene	42	ug/kg	1	-	U	Yes
1,2,4,5-Tetrachlorobenzene	210	ug/kg	1	-	U	Yes

METHOD: 8270D (SIM)

Naphthalene	4.2	ug/kg	1	-	U	Yes
1,4-Dioxane	8.07	ug/kg	1	-	-	Yes

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
--------------	--------	-------	-----------------	----------	------------	------------

Sample ID: JC19644-5

Sample location: BMSMC Building 5 Area

Sampling date: 5/4/2016

Matrix: Soil

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
2-Chlorophenol	81	ug/kg	1	-	U	Yes
4-Chloro-3-methyl phenol	200	ug/kg	1	-	U	Yes
2,4-Dichlorophenol	200	ug/kg	1	-	U	Yes
2,4-Dimethylphenol	200	ug/kg	1	-	U	Yes
2,4-Dinitrophenol	200	ug/kg	1	-	U	Yes
4,6-Dinitro-o-cresol	200	ug/kg	1	-	U	Yes
2-Methylphenol	81	ug/kg	1	-	U	Yes
3&4-Methylphenol	81	ug/kg	1	-	U	Yes
2-Nitrophenol	200	ug/kg	1	-	U	Yes
4-Nitrophenol	400	ug/kg	1	-	U	Yes
Pentachlorophenol	200	ug/kg	1	-	U	Yes
Phenol	81	ug/kg	1	-	U	Yes
2,3,4,6-Tetrachlorophenol	200	ug/kg	1	-	U	Yes
2,4,5-Trichlorophenol	200	ug/kg	1	-	U	Yes
2,4,6-Trichlorophenol	200	ug/kg	1	-	U	Yes
Acenaphthene	40	ug/kg	1	-	U	Yes
Acenaphthylene	40	ug/kg	1	-	U	Yes
Acetophenone	200	ug/kg	1	-	U	Yes
Anthracene	40	ug/kg	1	-	U	Yes
Atrazine	81	ug/kg	1	-	U	Yes
Benzo(a)anthracene	40	ug/kg	1	-	U	Yes
Benzo(a)pyrene	40	ug/kg	1	-	U	Yes
Benzo(b)fluoranthene	40	ug/kg	1	-	U	Yes
Benzo(g,h,i)perylene	40	ug/kg	1	-	U	Yes
Benzo(k)fluoranthene	40	ug/kg	1	-	U	Yes
4-Bromophenyl phenyl ether	81	ug/kg	1	-	U	Yes
Butyl benzyl phthalate	81	ug/kg	1	-	U	Yes
1,1'-Biphenyl	81	ug/kg	1	-	U	Yes
Benzaldehyde	200	ug/kg	1	-	U	Yes
2-Chloronaphthalene	81	ug/kg	1	-	U	Yes
4-Chloroaniline	200	ug/kg	1	-	U	Yes
Carbazole	81	ug/kg	1	-	U	Yes
Caprolactam	81	ug/kg	1	-	U	Yes
Chrysene	40	ug/kg	1	-	U	Yes
bis(2-Chloroethoxy)methane	81	ug/kg	1	-	U	Yes
bis(2-Chloroethyl)ether	81	ug/kg	1	-	U	Yes

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
bis(2-Chloroisopropyl)ether	81	ug/kg	1	-	U	Yes
4-Chlorophenyl phenyl ether	81	ug/kg	1	-	U	Yes
2,4-Dinitrotoluene	40	ug/kg	1	-	U	Yes
2,6-Dinitrotoluene	40	ug/kg	1	-	U	Yes
3,3'-Dichlorobenzidine	81	ug/kg	1	-	U	Yes
Dibenzo(a,h)anthracene	40	ug/kg	1	-	U	Yes
Dibenzofuran	81	ug/kg	1	-	U	Yes
Di-n-butyl phthalate	81	ug/kg	1	-	U	Yes
Di-n-octyl phthalate	81	ug/kg	1	-	U	Yes
Diethyl phthalate	81	ug/kg	1	-	U	Yes
Dimethyl phthalate	81	ug/kg	1	-	U	Yes
bis(2-Ethylhexyl)phthalate	81	ug/kg	1	-	U	Yes
Fluoranthene	40	ug/kg	1	-	U	Yes
Fluorene	40	ug/kg	1	-	U	Yes
Hexachlorobenzene	81	ug/kg	1	-	U	Yes
Hexachlorobutadiene	40	ug/kg	1	-	UJ	Yes
Hexachlorocyclopentadiene	400	ug/kg	1	-	U	Yes
Hexachloroethane	200	ug/kg	1	-	U	Yes
Indeno(1,2,3-cd)pyrene	40	ug/kg	1	-	U	Yes
Isophorone	81	ug/kg	1	-	U	Yes
1-Methylnaphthalene	81	ug/kg	1	-	U	Yes
2-Methylnaphthalene	81	ug/kg	1	-	U	Yes
2-Nitroaniline	200	ug/kg	1	-	U	Yes
3-Nitroaniline	200	ug/kg	1	-	U	Yes
4-Nitroaniline	200	ug/kg	1	-	U	Yes
Nitrobenzene	81	ug/kg	1	-	U	Yes
N-Nitroso-di-n-propylamine	81	ug/kg	1	-	U	Yes
Nitrosodiphenylamine	200	ug/kg	1	-	U	Yes
Phenanthrene	40	ug/kg	1	-	U	Yes
Pyrene	40	ug/kg	1	-	U	Yes
1,2,4,5-Tetrachlorobenzene	200	ug/kg	1	-	UJ	Yes

METHOD: 8270D (SIM)

Naphthalene	4.0	ug/kg	1	-	U	Yes
1,4-Dioxane	136	ug/kg	1	-	-	Yes

DATA REVIEW WORKSHEETS

Project Number: JC19644
 Date: May 03-04, 2016
 Shipping Date: May 03, 2016
 EPA Region: 2

REVIEW OF SEMIVOLATILE ORGANIC PACKAGE

The following guidelines for evaluating volatile organics were created to delineate required validation actions. This document will assist the reviewer in using professional judgment to make more informed decision and in better serving the needs of the data users. The sample results were assessed according to USEPA data validation guidance documents in the following order of precedence: EPA Hazardous Waste Support Section, SOP HW-35A, July 2015 –Revision 0. *Semivolatiles Data Validation*. The QC criteria and data validation actions listed on the data review worksheets are from the primary guidance document, unless otherwise noted.

The hardcopied (laboratory name) Accutest data package received has been reviewed and the quality control and performance data summarized. The data review for SVOCs included:

Lab. Project/SDG No.: JC19644 Sample matrix: Soil/Groundwater
 No. of Samples: 5 Full scan/5 SIM

Trip blank No.:
 Field blank No.:
 Equipment blank No.:
 Field duplicate No.:

<input checked="" type="checkbox"/> Data Completeness	<input checked="" type="checkbox"/> Laboratory Control Spikes
<input checked="" type="checkbox"/> Holding Times	<input checked="" type="checkbox"/> Field Duplicates
<input checked="" type="checkbox"/> GC/MS Tuning	<input checked="" type="checkbox"/> Calibrations
<input checked="" type="checkbox"/> Internal Standard Performance	<input checked="" type="checkbox"/> Compound Identifications
<input checked="" type="checkbox"/> Blanks	<input checked="" type="checkbox"/> Compound Quantitation
<input checked="" type="checkbox"/> Surrogate Recoveries	<input checked="" type="checkbox"/> Quantitation Limits
<input checked="" type="checkbox"/> Matrix Spike/Matrix Spike Duplicate	

Overall Comments: ABN TCL list by method SW846-8270D; Naphthalene and 1,4-Dioxane analyzed by method SW846-8270D (SIM)

Definition of Qualifiers:

J- Estimated results
 U- Compound not detected
 R- Rejected data
 UJ- Estimated nondetect

Reviewer: Rafael Infante
 Date: May 30, 2016

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

HOLDING TIMES

The objective of this parameter is to ascertain the validity of the results based on the holding time of the sample from time of collection to the time of analysis.

Complete table for all samples and note the analysis and/or preservation not within criteria

SAMPLE ID	DATE SAMPLED	DATE EXTRACTED/ANALYZED	pH	ACTION
All samples extracted and analyzed within method recommended holding time. Sample preservation was acceptable.				

Cooler temperature (Criteria: 4 ± 2 °C): 4°C

Actions

Results will be qualified based on the criteria of the following Table:

Table 1. Holding Time Actions for Semivolatile Analyses

Matrix	Preserved	Criteria	Action	
			Detected Associated Compounds	Non-Detected Associated Compounds
Aqueous	No	≤ 7 days (for extraction) ≤ 40 days (for analysis)	Use professional judgment	
	No	> 7 days (for extraction) > 40 days (for analysis)	J	Use professional judgment
	Yes	≤ 7 days (for extraction) ≤ 40 days (for analysis)	No qualification	
	Yes	> 7 days (for extraction) > 40 days (for analysis)	J	UJ
	Yes/No	Grossly Exceeded	J	UJ or R
Non-Aqueous	No	≤ 14 days (for extraction) ≤ 40 days (for analysis)	Use professional judgment	
	No	> 14 days (for extraction) > 40 days (for analysis)	J	Use professional judgment
	Yes	≤ 14 days (for extraction) ≤ 40 days (for analysis)	No qualification	
	Yes	> 14 days (for extraction) > 40 days (for analysis)	J	UJ
	Yes/No	Grossly Exceeded	J	UJ or R

DATA REVIEW WORKSHEETS

All criteria were met X
Criteria were not met see below

GC/MS TUNING

The assessment of the tuning results is to determine if the sample instrumentation is within the standard tuning QC limits

 X The DFTPP performance results were reviewed and found to be within the specified criteria.

 X DFTPP tuning was performed for every 12 hours of sample analysis.

If no, use professional judgment to determine whether the associated data should be accepted, qualified or rejected.

Notes: These requirements do not apply when samples are analyzed by the Selected Ion Monitoring (SIM) technique.

All mass spectrometer conditions must be identical to those used during the sample analysis. Background subtraction actions resulting in spectral distortion are unacceptable

Notes: No data should be qualified based of DFTPP failure.

The requirement to analyze the instrument performance check solution is optional when analysis of PAHs/pentachlorophenol is to be performed by the SIM technique.

List	the	samples	affected:
_____			_____
_____			_____
_____			_____
_____			_____

Actions:

1. If sample are analyzed without a preceding valid instrument performance check or are analyzed 12 hours after the Instrument Performance Check, qualify all data in those samples as unusable (R).
2. If ion abundance criteria are not met, use professional judgment to determine to what extent the data may be utilized.
3. State in the Data Review Narrative, decisions to use analytical data associated with DFTPP instrument performance checks not meeting the contract requirements.
4. Use professional judgment to determine if associated data should be qualified based on the spectrum of the mass calibration compounds.

DATA REVIEW WORKSHEETS

All criteria were met ☒ X
 Criteria were not met
 and/or see below _____

INITIAL CALIBRATION VERIFICATION

Compliance requirements for satisfactory instrument calibration are established to ensure that the instrument is capable of producing and maintaining acceptable quantitative data.

Date of initial calibration: _____ 04/21/2016 (SIM) _____ 04/04-05/2016 (Scan) _____
 Instrument ID numbers: _____ GCMS3M _____ GCMSF _____
 Matrix/Level: _____ Aqueous/low _____ Aqueous/low _____
 Date of initial calibration: _____ 04/26/16; 04/27/16 (Scan) _____ 04/14/16 (SIM) _____
 Instrument ID numbers: _____ GCMS3E _____ GCMS4M _____
 Matrix/Level: _____ Aqueous/low _____ Aqueous/low _____

DATE	LAB FILE ID#	CRITERIA OUT RFs, %RSD, %D, r	COMPOUND	SAMPLES AFFECTED
Initial and initial calibration verification meets the method and guidance validation document performance criteria.				

Note: GCMS6P instrument used in the scan mode. Calibrations were within guidance document criteria. Used for the analysis of QC sample only.

Actions:

Qualify the initial calibration analytes listed in Table 2 using the following criteria:

Table 3. Initial Calibration Actions for Semivolatile Analysis

Criteria	Action	
	Detect	Non-detect
Initial Calibration not performed at specified frequency and sequence	Use professional judgment R	Use professional judgment R
Initial Calibration not performed at the specified concentrations	J	UJ
RRF < Minimum RRF in Table 2 for target analyte	Use professional judgment J+ or R	R
RRF ≥ Minimum RRF in Table 2 for target analyte	No qualification	No qualification
%RSD > Maximum %RSD in Table 2 for target analyte	J	Use professional judgment
%RSD ≤ Maximum %RSD in Table 2 for target analyte	No qualification	No qualification

DATA REVIEW WORKSHEETS

Initial Calibration

Table 2. RRF, %RSD, and %D Acceptance Criteria in Initial Calibration and CCV for Semivolatili Analysis

Analyte	Minimum RRF	Maximum %RSD	Opening Maximum %D ¹	Opening Maximum %D ¹
1,4-Dioxane	0.010	40.0	± 40.0	± 50.0
Benzaldehyde	0.100	40.0	± 40.0	± 50.0
Phenol	0.080	20.0	± 20.0	± 25.0
Bis(2-chloroethyl)ether	0.100	20.0	± 20.0	± 25.0
2-Chlorophenol	0.200	20.0	± 20.0	± 25.0
2-Methylphenol	0.010	20.0	± 20.0	± 25.0
3-Methylphenol	0.010	20.0	± 20.0	± 25.0
2,2'-Oxybis-(1-chloropropane)	0.010	20.0	± 25.0	± 50.0
Acetophenone	0.060	20.0	± 20.0	± 25.0
4-Methylphenol	0.010	20.0	± 20.0	± 25.0
N-Nitroso-di-n-propylamine	0.080	20.0	± 25.0	± 25.0
Hexachloroethane	0.100	20.0	± 20.0	± 25.0
Nitrobenzene	0.090	20.0	± 20.0	± 25.0
Isophorone	0.100	20.0	± 20.0	± 25.0
2-Nitrophenol	0.060	20.0	± 20.0	± 25.0
2,4-Dimethylphenol	0.050	20.0	± 25.0	± 50.0
Bis(2-chloroethoxy)methane	0.080	20.0	± 20.0	± 25.0
2,4-Dichlorophenol	0.060	20.0	± 20.0	± 25.0
Naphthalene	0.200	20.0	± 20.0	± 25.0
4-Chloroaniline	0.010	40.0	± 40.0	± 50.0
Hexachlorobutadiene	0.040	20.0	± 20.0	± 25.0
Caprolactam	0.010	40.0	± 30.0	± 50.0
4-Chloro-3-methylphenol	0.040	20.0	± 20.0	± 25.0
2-Methylnaphthalene	0.100	20.0	± 20.0	± 25.0
Hexachlorocyclopentadiene	0.010	40.0	± 40.0	± 50.0
2,4,6-Trichlorophenol	0.090	20.0	± 20.0	± 25.0
2,4,5-Trichlorophenol	0.100	20.0	± 20.0	± 25.0
1,1'-Biphenyl	0.200	20.0	± 20.0	± 25.0

DATA REVIEW WORKSHEETS

Analyte	Minimum RRF	Maximum %RSD	Opening Maximum %D ¹	Opening Maximum %D ¹
2-Chloronaphthalene	0.300	20.0	± 20.0	± 25.0
2-Nitroaniline	0.060	20.0	± 25.0	± 25.0
Dimethylphthalate	0.300	20.0	± 25.0	± 25.0
2,6-Dinitrotoluene	0.080	20.0	± 20.0	± 25.0
Acenaphthylene	0.400	20.0	± 20.0	± 25.0
3-Nitroaniline	0.010	20.0	± 25.0	± 50.0
Acenaphthene	0.200	20.0	± 20.0	± 25.0
2,4-Dinitrophenol	0.010	40.0	± 50.0	± 50.0
4-Nitrophenol	0.010	40.0	± 40.0	± 50.0
Dibenzofuran	0.300	20.0	± 20.0	± 25.0
2,4-Dinitrotoluene	0.070	20.0	± 20.0	± 25.0
Diethylphthalate	0.300	20.0	± 20.0	± 25.0
1,2,4,5-Tetrachlorobenzene	0.100	20.0	± 20.0	± 25.0
4-Chlorophenyl-phenylether	0.100	20.0	± 20.0	± 25.0
Fluorene	0.200	20.0	± 20.0	± 25.0
4-Nitroaniline	0.010	40.0	± 40.0	± 50.0
4,6-Dinitro-2-methylphenol	0.010	40.0	± 30.0	± 50.0
4-Bromophenyl-phenyl ether	0.070	20.0	± 20.0	± 25.0
N-Nitrosodiphenylamine	0.100	20.0	± 20.0	± 25.0
Hexachlorobenzene	0.050	20.0	± 20.0	± 25.0
Atrazine	0.010	40.0	± 25.0	± 50.0
Pentachlorophenol	0.010	40.0	± 40.0	± 50.0
Phenanthrene	0.200	20.0	± 20.0	± 25.0
Anthracene	0.200	20.0	± 20.0	± 25.0
Carbazole	0.050	20.0	± 20.0	± 25.0
Di-n-butylphthalate	0.500	20.0	± 20.0	± 25.0
Fluoranthene	0.100	20.0	± 20.0	± 25.0
Pyrene	0.400	20.0	± 25.0	± 50.0
Butylbenzylphthalate	0.100	20.0	± 25.0	± 50.0

DATA REVIEW WORKSHEETS

Analyte	Minimum RRF	Maximum %RSD	Opening Maximum %D'	Opening Maximum %D'
3,3'-Dichlorobenzidine	0.010	40.0	± 40.0	± 50.0
Benzo(a)anthracene	0.300	20.0	± 20.0	± 25.0
Chrysene	0.200	20.0	± 20.0	± 50.0
Bis(2-ethylhexyl) phthalate	0.200	20.0	± 25.0	± 50.0
Di-n-octylphthalate	0.010	40.0	± 40.0	± 50.0
Benzo(b)fluoranthene	0.010	20.0	± 25.0	± 50.0
Benzo(k)fluoranthene	0.010	20.0	± 25.0	± 50.0
Benzo(a)pyrene	0.010	20.0	± 20.0	± 50.0
Indeno(1,2,3-cd)pyrene	0.010	20.0	± 25.0	± 50.0
Dibenzo(a,h)anthracene	0.010	20.0	± 25.0	± 50.0
Benzo(g,h,i)perylene	0.010	20.0	± 30.0	± 50.0
2,3,4,6-Tetrachlorophenol	0.040	20.0	± 20.0	± 50.0
Naphthalene	0.600	20.0	± 25.0	± 25.0
2-Methylnaphthalene	0.300	20.0	± 20.0	± 25.0
Acenaphthylene	0.900	20.0	± 20.0	± 25.0
Acenaphthene	0.500	20.0	± 20.0	± 25.0
Fluorene	0.700	20.0	± 25.0	± 50.0
Phenanthrene	0.300	20.0	± 25.0	± 50.0
Anthracene	0.400	20.0	± 25.0	± 50.0
Fluoranthene	0.400	20.0	± 25.0	± 50.0
Pyrene	0.500	20.0	± 30.0	± 50.0
Benzo(a)anthracene	0.400	20.0	± 25.0	± 50.0
Chrysene	0.400	20.0	± 25.0	± 50.0
Benzo(b)fluoranthene	0.100	20.0	± 30.0	± 50.0
Benzo(k)fluoranthene	0.100	20.0	± 30.0	± 50.0
Benzo(a)pyrene	0.100	20.0	± 25.0	± 50.0
Indeno(1,2,3-cd)pyrene	0.100	20.0	± 40.0	± 50.0
Dibenzo(a,h)anthracene	0.010	25.0	± 40.0	± 50.0
Benzo(g,h,i)perylene	0.020	25.0	± 40.0	± 50.0

DATA REVIEW WORKSHEETS

Pentachlorophenol	0.010	40.0	± 50.0	± 50.0
Deuterated Monitoring Compounds				
Analyte	Minimum RRF	Maximum %RSD	Opening Maximum %D ¹	Closing Maximum %D
1,4-Dioxane-d ₈	0.010	20.0	± 25.0	± 50.0
Phenol-d ₅	0.010	20.0	± 25.0	± 25.0
Bis-(2-chloroethyl)ether-d ₈	0.100	20.0	± 20.0	± 25.0
2-Chlorophenol-d ₄	0.200	20.0	± 20.0	± 25.0
4-Methylphenol-d ₈	0.010	20.0	± 20.0	± 25.0
4-Chloroaniline-d ₄	0.010	40.0	± 40.0	± 50.0
Nitrobenzene-d ₅	0.050	20.0	± 20.0	± 25.0
2-Nitrophenol-d ₄	0.050	20.0	± 20.0	± 25.0
2,4-Dichlorophenol-d ₃	0.060	20.0	± 20.0	± 25.0
Dimethylphthalate-d ₆	0.300	20.0	± 20.0	± 25.0
Acenaphthylene-d ₈	0.400	20.0	± 20.0	± 25.0
4-Nitrophenol-d ₄	0.010	40.0	± 40.0	± 50.0
Fluorene-d ₁₀	0.100	20.0	± 20.0	± 25.0
4,6-Dinitro-2-methylphenol-d ₂	0.010	40.0	± 30.0	± 50.0
Anthracene-d ₁₀	0.300	20.0	± 20.0	± 25.0
Pyrene-d ₁₀	0.300	20.0	± 25.0	± 50.0
Benzo(a)pyrene-d ₁₂	0.010	20.0	± 20.0	± 50.0
Fluoranthene-d ₁₀ (SIM)	0.400	20.0	± 25.0	± 50.0
2-Methylnaphthalene-d ₁₀ (SIM)	0.300	20.0	± 20.0	± 25.0

¹ If a closing CCV is acting as an opening CCV, all target analytes must meet the requirements for an opening CCV.

Note: If analysis by SIM technique is requested for PAH/pentachlorophenols, calibration standards analyzed at 0.10, 0.20, 0.40, 0.80, and 1.0 ng/uL for each target compound of interest and the associated DMCs. Pentachlorophenol will require only a four point initial calibration at 0.20, 0.40, 0.80, and 1.0 ng/uL.

DATA REVIEW WORKSHEETS

All criteria were met _____
 Criteria were not met _____
 and/or see below ____X____

CONTINUING CALIBRATION VERIFICATION

Compliance requirements for satisfactory instrument calibration are established to ensure that the instrument is capable of producing and maintaining acceptable quantitative data.

Date of initial calibration: _____ 04/04-05/16 _____ (Scan) _____
 Date of initial calibration verification (ICV): _____ 04/05-06/16 _____
 Date of continuing calibration verification (CCV): _____ 05/06/16 _____
 Date of closing CCV: _____ - _____
 Instrument ID numbers: _____ GCMSF _____
 Matrix/Level: _____ Aqueous/low _____

Date of initial calibration: _____ 04/21/16 _____ (SIM) _____
 Date of initial calibration verification (ICV): _____ 04/21/16 _____
 Date of continuing calibration verification (CCV): _____ 05/06/16 _____
 Date of closing CCV: _____ - _____
 Instrument ID numbers: _____ GCMS3M _____
 Matrix/Level: _____ Aqueous/low _____

Date of initial calibration: _____ 04/26/16; 04/27/16 _____ (Scan) _____
 Date of initial calibration verification (ICV): _____ 04/26-27/16; 04/27/16 _____
 Date of continuing calibration verification (CCV): _____ 05/06/16; 05/09/11 _____
 Date of closing CCV: _____ - _____
 Instrument ID numbers: _____ GCMS3E _____
 Matrix/Level: _____ Aqueous/low _____

Date of initial calibration: _____ 04/14/16 _____ (SIM) _____
 Date of initial calibration verification (ICV): _____ 04/14/16 _____
 Date of continuing calibration verification (CCV): _____ 05/09/16 _____
 Date of closing CCV: _____ - _____
 Instrument ID numbers: _____ GCMS4M _____
 Matrix/Level: _____ Aqueous/low _____

DATE	LAB FILE ID#	CRITERIA OUT RFs, %RSD, %D, r	COMPOUND	SAMPLES AFFECTED
GCMS3E				
05/06/16	cc3623-50	-21.6	Hexachlorobutadiene	JC19644-4; -5
		-22.8	Pentachlorophenol*	
05/06/16	cc3621-50	32.5	Benzaldehyde*	
		-24.8	1,2,4,5-tetrachlorobenzene	
GCMSF				
05/06/16	cc6563-25	-26.9	n-Nitroso-di-n-propylamine	JC19644-1; -2; -3
		-25.8	2-Nitroaniline	

DATA REVIEW WORKSHEETS

Note: Initial and continuing calibration verifications meet the method and guidance document required performance criteria except the cases described above. Results qualified as estimated (J), (UJ) for non-detects.

No closing calibration verification included in data package. No action taken, professional judgment.

* Analytes with % difference in the continue calibration verification outside the method performance criteria but within the validation guidelines criteria, + 40 %. No action taken.

QC samples for this data package were analyzed in the instrument GCM6P (Scan). Analytes not meeting continue calibration verification criteria were not qualified.

Actions:

Notes: Verify that the CCV is run at the required frequency (an opening and closing CCV must be run within 12-hour period).

All DMCs must meet the RRF values given in Table 2. No qualification of the data is necessary on DMCs RRF and %RSD/%D alone. Use professional judgment to evaluate DMCs and %RSD/%D data in conjunction with DMCs recoveries to determine the need for qualification of the data.

Qualify the initial calibration analytes listed in Table 2 using the following criteria in the CCVs:

Table 4. CCV Actions for Semivolatile Analysis

Criteria for Opening CCV	Criteria for Closing CCV	Action	
		Detect	Non-detect
CCV not performed at required frequency and sequence	CCV not performed at required frequency	Use professional judgment R	Use professional judgment R
CCV not performed at specified concentration	CCV not performed at specified concentration	Use professional judgment	Use professional judgment
RRF < Minimum RRF in Table 2 for target analyte	RRF < Minimum RRF in Table 2 for target analyte	Use professional judgment J or R	R
RRF ≥ Minimum RRF in Table 2 for target analyte	RRF ≥ Minimum RRF in Table 2 for target analyte	No qualification	No qualification
%D outside the Opening Maximum %D limits in Table 2 for target analyte	%D outside the Closing Maximum %D limits in Table 2 for target analyte	J	UJ
%D within the inclusive Opening Maximum %D limits in Table 2 for target analyte	%D within the inclusive Closing Maximum %D limits in Table 2 for target analyte	No qualification	No qualification

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

BLANK ANALYSIS RESULTS (Sections 1 & 2)

The assessment of the blank analysis results is to determine the existence and magnitude of contamination problems. The criteria for evaluation of blanks apply only to blanks associated with the samples, including trip, equipment, and laboratory blanks. If problems with any blanks exist, all data associated with the case must be carefully evaluated to determine whether or not there is an inherent variability in the data for the case, or if the problem is an isolated occurrence not affecting other data.

List the contamination in the blanks below. High and low levels blanks must be treated separately.

Notes: The concentration of non-target compounds in all blanks must be less than or equal to 10 ug/L.

The concentration of target compounds in all blanks must be less than its CRQL listed in the method.

Samples taken from a drinking water tap do not have an associated field blank.

Laboratory blanks

DATE ANALYZED	LAB ID	LEVEL/ MATRIX	COMPOUND	CONCENTRATION UNITS
_No_target_analytes_detected_in_method_blanks._				

Field/Equipment/Trip blank

DATE ANALYZED	LAB ID	LEVEL/ MATRIX	COMPOUND	CONCENTRATION UNITS
_No_field/trip/equipment_blanks_analyzed_with_this_data_package._				

DATA REVIEW WORKSHEETS

All criteria were met ☒ X
 Criteria were not met
 and/or see below _____

BLANK ANALYSIS RESULTS (Section 3)

Blank Actions

Qualify samples based on the criteria summarized in Table 5:

Table 5. Blank and TCLP/SPLP LEB Actions for Semivolatile Analysis

Blank Type	Blank Result	Sample Result	Action
Method, TCLP/SPLP LEB, Field	Detect	Non-detect	No qualification
	< CRQL	< CRQL	Report at CRQL and qualify as non-detect (U)
		≥ CRQL	Use professional judgment
	≥ CRQL	< CRQL	Report at CRQL and qualify as non-detect (U)
		≥ CRQL but < Blank Result	Report at sample results and qualify as non-detect (U) or as unusable (R)
		≥ CRQL and ≥ Blank Result	Use professional judgment
	Grossly high	Detect	Report at sample results and qualify as unusable (R)
	TIC > 5.0 ug/L (water) or 0.0050 mg/L (TCLP leachate) or TIC > 170 ug/Kg (soil)	Detect	Use professional judgment

List samples qualified

CONTAMINATION SOURCE/LEVEL	COMPOUND	CONC/UNITS	AL/UNITS	SQL	AFFECTED SAMPLES

All criteria were met _____
 Criteria were not met _____
 and/or see below X

SURROGATE SPIKE RECOVERIES – DEUTERATED MONITORING COMPOUNDS (DMCs)

Laboratory performance of individual samples is established by evaluation of surrogate spike recoveries – deuterated monitoring compounds. All samples are spiked with surrogate compounds prior to sample analysis. The accuracy of the analysis is measured by the surrogate percent recovery. Since the effects of the sample matrix are frequently outside the control of the laboratory and may present relatively unique problems, the validation of data is frequently subjective and demands analytical experience and professional judgment.

Notes: Recoveries for DMCs in samples and blanks must be within the limits specified in Table 6.

The recovery limits for any of the compounds listed in Table 6 may be expanded at any time during the period of performance if USEPA determines that the limits are too restrictive.

If a DMC is not added in the samples and blanks or the concentrations of DMCs in the samples and blank not the specified, use professional judgment in qualifying the data.

Table 7. DMC Actions for Semivolatile Analysis

Criteria	Action	
	Detect	Non-detect
%R < 10% (excluding DMCs with 10% as a lower acceptance limit)	J-	R
10% ≤ %R (excluding DMCs with 10% as a lower acceptance limit) < Lower Acceptance Limit	J-	UJ
Lower Acceptance limit ≤ %R ≤ Upper Acceptance Limit	No qualification	No qualification
%R > Upper Acceptance Limit	J+	No qualification

List the percent recoveries (%Rs) which do not meet the criteria for DMCs (surrogate) recovery.

Matrix: Groundwater

SAMPLE ID	SURROGATE COMPOUND	ACTION
-----------	--------------------	--------

DMCs meet the required criteria except for the cases described below. Non-deuterated
surrogates added to the samples were within laboratory recovery limits unless the cases
described in this document.

Samples and QC shown here apply to the above method

DATA REVIEW WORKSHEETS

Lab Sample ID	Lab File ID	S1	S2	S3	S4	S5	S6
JC19644-2	F156880.D	0* a	0* a	0* a	70	63	71
Surrogate Compounds	Recovery Limits	Surrogate Compounds	Recovery Limits				
S1 = 2-Fluorophenol	14-88%	S4 = Nitrobenzene-d5	32-128%				
S2 = Phenol-d5	10-110%	S5 = 2-Fluorobiphenyl	35-119%				
S3 = 2,4,6-Tribromophenol	39-149%	S6 = Terphenyl-d14	10-126%				

(a) Outside control limits due to dilution. No action taken.

Table 8. Semivolatile DMCs and the Associated Target Analytes

1,4-Dioxane-d ₈ (DMC-1)	Phenol-d ₅ (DMC-2)	Bis(2-Chloroethyl) ether-d ₈ (DMC-3)
1,4-Dioxane	Benzaldehyde Phenol	Bis(2-chloroethyl) ether 2,2'-Oxybis(1-chloropropane) Bis(2-chloroethoxy)methane
2-Chlorophenol-d ₄ (DMC-4)	4-Methylphenol-d ₄ (DMC-5)	4-Chloroaniline-d ₄ (DMC-6)
2-Chlorophenol	2-Methylphenol 3-Methylphenol 4-Methylphenol 2,4-Dimethylphenol	4-Chloroaniline Hexachlorocyclopentadiene Dichlorobenzidine
Nitrobenzene-d ₅ (DMC-7)	2-Nitrophenol-d ₄ (DMC-8)	2,4-Dichlorophenol-d ₃ (DMC-9)
Acetophenone N-Nitroso-di-n-propylamine Hexachloroethane Nitrobenzene 2,6-Dinitrotoluene 2,4-Dinitrotoluene N-Nitrosodiphenylamine	Isophorone 2-Nitrophenol	2,4-Dichlorophenol Hexachlorobutadiene Hexachlorocyclopentadiene 4-Chloro-3-methylphenol 2,4,6-Trichlorophenol 2,4,5-Trichlorophenol 1,2,4,5-Tetrachlorobenzene *Pentachlorophenol 2,3,4,6-Tetrachlorophenol
Dimethylphthalate-d ₄ (DMC-10)	Acenaphthylene-d ₈ (DMC-11)	4-Nitrophenol-d ₄ (DMC-12)
Caprolactam 1,1'-Biphenyl Dimethylphthalate Diethylphthalate Di-n-butylphthalate Butylbenzylphthalate Bis(2-ethylhexyl) phthalate Di-n-octylphthalate	*Naphthalene *2-Methylnaphthalene 2-Chloronaphthalene *Acenaphthylene *Acenaphthene	2-Nitroaniline 3-Nitroaniline 2,4-Dinitrophenol 4-Nitrophenol 4-Nitroaniline

DATA REVIEW WORKSHEETS

Fluorene-d ₁₀ (DMC-13)	4,6-Dinitro-2-methylphenol-d ₂ (DMC-14)	Anthracene-d ₁₀ (DMC-15)
Dibenzofuran *Fluorene 4-Chlorophenyl-phenylether 4-Bromophenyl-phenylether Carbazole	4,6-Dinitro-2-methylphenol	Hexachlorobenzene Atrazine *Phenanthrene *Anthracene
Pyrene-d ₁₀ (DMC-16)	Benzo(a)pyrene-d ₁₂ (DMC-17)	
*Fluoranthene *Pyrene *Benzo(a)anthracene *Chrysene	3,3'-Dichlorobenzidine *Benzo(b)fluoranthene *Benzo(k)fluoranthene *Benzo(a)pyrene *Indeno(1,2,3-cd)pyrene *Dibenzo(a,h)anthracene *Benzo(g,h,i)perylene	

*Included in optional Target Analyte List (TAL) of PAHs and PCP only.

Table 9. Semivolatile SIM DMCs and the Associated Target Analytes

Fluoranthene-d ₁₀ (DMC-1)	2-Methylnaphthalene-d ₁₀ (DMC-2)
Fluoranthene	Naphthalene
Pyrene	2-Methylnaphthalene
Benzo(a)anthracene	Acenaphthylene
Chrysene	Acenaphthene
Benzo(b)fluoranthene	Fluorene
Benzo(k)fluoranthene	Pentachlorophenol
Benzo(a)pyrene	Phenanthrene
Indeno(1,2,3-cd)pyrene	Anthracene
Dibenzo(a,h)anthracene	
Benzo(g,h,i)perylene	

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met _____
 and/or see below _____

VII. A MATRIX SPIKE/MATRIX SPIKE DUPLICATE (MS/MSD)

This data is generated to determine long term precision and accuracy in the analytical method for various matrices. This data alone cannot be used to evaluate the precision and accuracy of individual samples. If any % R in the MS or MSD falls outside the designated range, the reviewer should determine if there are matrix effects, i.e. LCS data are within the QC limits but MS/MSD data are outside QC limit.

1. MS/MSD Recoveries and Precision Criteria

The laboratory should use one MS and a duplicate analysis of an unspiked field sample if target analytes are expected in the sample. If target analytes are not expected, MS/MSD should be analyzed.

NOTES: Data for MS and MSDs will not be present unless requested by the Region.
 Notify the Contract Laboratory COR if a field or trip blank was used for the MS and MSD.

For a Matrix Spike that does not meet criteria, apply the action to only the field sample used to prepare the Matrix Spike sample. If it is clearly stated in the data validation materials that the samples were taken through incremental sampling or some other method guaranteeing the homogeneity of the sample group, then the entire sample group may be qualified.

List the %Rs, RPD of the compounds which do not meet the criteria.

Sample ID: <u> JC19233-1 </u>	Matrix/Level: <u> Groundwater </u>
Sample ID: <u> JC19639-1 (SIM) </u>	Matrix/Level: <u> Groundwater </u>
Sample ID: <u> JC19596-12A </u>	Matrix/Level: <u> Soil </u>
Sample ID: <u> JC19644-4 (SIM) </u>	Matrix/Level: <u> Soil </u>

* see enclosed list

MS OR MSD	COMPOUND	% R	RPD	QC LIMITS	ACTION
<u> JC19596-12AMS/MSD </u>					
<u> MS/MSD </u>	<u> Hexachlorocyclopentadiene </u>	<u> 0/0 </u>	<u> % </u>	<u> 10 - 127 </u>	<u> No action </u>

Note: No action taken, MS/MSD results apply only to unspiked sample.

- * QC limits are laboratory in-house performance criteria, LL = lower limit, UL = upper limit.
- * If QC limits are not available, use limits of 70 – 130 %.

DATA REVIEW WORKSHEETS

Actions:

QUALITY	%R < LL	%R > UL
Positive results	J	J
Nondetects results	R	Accept

MS/MSD criteria apply only to the unspiked sample, its dilutions, and the associated MS/MSD samples:

If the % R for the affected compounds were < LL (or 70 %), qualify positive results (J) and nondetects (UJ).

If the % R for the affected compounds were > UL (or 130 %), only qualify positive results (J).

If 25 % or more of all MS/MSD %R were < LL (or 70 %) or if two or more MS/MSD %Rs were < 10%, qualify all positive results (J) and reject nondetects (R).

A separate worksheet should be used for each MS/MSD pair.

DATA REVIEW WORKSHEETS

All criteria were met ☒ X
 Criteria were not met
 and/or see below _____

INTERNAL STANDARD PERFORMANCE

The assessment of the internal standard (IS) parameter is used to assist the data reviewer in determining the condition of the analytical instrumentation.

List the internal standard area of samples which do not meet the criteria.

DATE	SAMPLE ID	IS OUT	IS AREA	ACCEPTABLE RANGE	ACTION
------	-----------	--------	---------	------------------	--------

Internal standard area counts meet the required criteria except for the cases described below.

05/19/16	JC19644-5	Perylene-d12	295560	300130	No action
05/19/16	JC19644-5*	Chrysene-d12	351314	363916	No action
05/19/16	JC19644-5*	Perylene-d12	290148	300130	No action

* Confirmation run for internal standard areas.

Action:

1. If an internal standard area count for a sample or blank is greater than 200.0% of the area for the associated standard (opening CCV or mid-point standard from initial calibration) (see Table 10 below):
 - a. Qualify detects for compounds quantitated using that internal standard as estimated low (J-).
 - b. Do not qualify non-detected associated compounds.
2. If an internal standard area count for a sample or blank is less than 20.0% of the area for the associated standard (opening CCV or mid-point standard from initial calibration):
 - a. Qualify detects for compounds quantitated using that internal standard as estimated high (J+).
 - b. Qualify non-detected associated compounds as unusable (R).
3. If an internal standard area count for a sample or blank is greater than or equal to 50.0%, and less than or equal to 200% of the area for the associated standard opening CCV or mid-point standard from initial calibration, no qualification of the data is necessary.
4. If an internal standard RT varies by more than 10.0 seconds: Examine the chromatographic profile for that sample to determine if any false positives or negatives exist. For shifts of a large magnitude, the reviewer may consider partial or total rejection of the data for that sample fraction. Detects should not need to be qualified as unusable (R) if the mass spectral criteria are met.
5. If an internal standard RT varies by less than or equal to 10.0 seconds, no qualification of the data is necessary.

Note: Inform the Contract Laboratory Program Project Officer (CLP PO) if the internal standard performance criteria are grossly exceeded. Note in the Data Review Narrative potential effects on the data resulting from unacceptable internal standard performance.

DATA REVIEW WORKSHEETS

State in the Data Review Narrative if the required internal standard compounds are not added to a sample or blank or if the required internal standard compound is not analyzed at the specified concentration.

Actions:

Table 10. Internal Standard Actions for Semivolatile Analysis

Criteria	Action	
	Detect	Non-detect
Area response < 20% of the opening CCV or mid-point standard CS3 from ICAL	J+	R
$20\% \leq$ Area response < 50% of the opening CCV or mid-point standard CS3 from ICAL	J+	UJ
$50\% \leq$ Area response \leq 200% of the opening CCV or mid-point standard CS3 from ICAL	No qualification	No qualification
Area response > 200% of the opening CCV or mid-point standard CS3 from ICAL	J-	No qualification
RT shift between sample/blank and opening CCV or mid-point standard CS3 from ICAL > 10.0 seconds	R	R
RT shift between sample/blank and opening CCV or mid-point standard CS3 from ICAL < 10.0 seconds	No qualification	No qualification

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

TARGET COMPOUND IDENTIFICATION

Criteria:

Is the Relative Retention Times (RRTs) of reported compounds within ± 0.06 RRT units of the standard RRT [opening Continuing Calibration Verification (CCV) or mid-point standard from the initial calibration]. **Yes? or No?**

List compounds not meeting the criteria described above:

Sample ID	Compounds	Actions
=====	=====	=====
_____	_____	_____
_____	_____	_____
_____	_____	_____

Mass spectra of the sample compound and a current laboratory-generated standard [i.e., the mass spectrum from the associated calibration standard (opening CCV or mid-point standard from initial calibration)] must match according to the following criteria:

- All ions present in the standard mass spectrum at a relative intensity greater than 10% must be present in the sample spectrum.
- The relative intensities of these ions must agree within $\pm 20\%$ between the standard and sample spectra (e.g., for an ion with an abundance of 50% in the standard spectrum, the corresponding sample ion abundance must be between 30-70%).
- Ions present at greater than 10% in the sample mass spectrum, but not present in the standard spectrum, must be evaluated by a reviewer experienced in mass spectral interpretation.

List compounds not meeting the criteria described above:

Sample ID	Compounds	Actions
=====	=====	=====
_____	_____	_____
_____	_____	_____
_____	_____	_____

Identified compounds meet the required criteria

DATA REVIEW WORKSHEETS

Action:

1. The application of qualitative criteria for GC/MS analysis of target compounds requires professional judgment. It is up to the reviewer's discretion to obtain additional information from the laboratory. If it is determined that incorrect identifications were made, qualify all such data as unusable (R).
2. Use professional judgment to qualify the data if it is determined that cross-contamination has occurred.
3. Note in the Data Review Narrative any changes made to the reported compounds or concerns regarding target compound identifications. Note, for Contract Laboratory COR action, the necessity for numerous or significant changes.

TENTATIVELY IDENTIFIED COMPOUNDS (TICS)

NOTE: Tentatively identified compounds should only be evaluated when requested by a party from outside of the Hazardous Waste Support Section (HWSS).

List TICs

Sample ID	Compound	Sample ID	Compound
=====			

Action:

1. Qualify all TIC results for which there is presumptive evidence of a match (e.g. greater than or equal to 85% match) as tentatively identified (NJ), with approximated concentrations. TICs labeled "unknown" are qualified as estimated (J).
2. General actions related to the review of TIC results are as follows:
 - a. If it is determined that a tentative identification of a non-target compound is unacceptable, change the tentative identification to "unknown" or another appropriate identification, and qualify the result as estimated (J).
 - b. If all contractually-required peaks were not library searched and quantitated, the Region's designated representative may request these data from the laboratory.
3. In deciding whether a library search result for a TIC represents a reasonable identification, use professional judgment. If there is more than one possible match, report the result as "either compound X or compound Y". If there is a lack of isomer specificity, change the TIC result to a nonspecific isomer result (e.g., 1,3,5-trimethyl benzene to trimethyl benzene isomer) or to a compound class (e.g., 2-methyl, 3-ethyl benzene to a substituted aromatic compound).
4. The reviewer may elect to report all similar compounds as a total (e.g., all alkanes may be summarized and reported as total hydrocarbons).

DATA REVIEW WORKSHEETS

5. Target compounds from other fractions and suspected laboratory contaminants should be marked as "non-reportable".
6. Other Case factors may influence TIC judgments. If a sample TIC match is poor, but other samples have a TIC with a valid library match, similar RRT, and the same ions, infer identification information from the other sample TIC results.
7. Note in the Data Review Narrative any changes made to the reported data or any concerns regarding TIC identifications.
8. Note, for Contract Laboratory COR action, failure to properly evaluate and report TICs

All criteria were met X
 Criteria were not met
 and/or see below

SAMPLE QUANTITATION AND REPORTED CONTRACT REQUIRED QUANTITATION LIMITS (CRQLS)

Action:

1. When a sample is analyzed at more than one dilution, the lower CRQL are used unless a QC exceedance dictates the use of higher CRQLs from the diluted sample. Samples reported with an "E" qualifier should be reported from the diluted sample.
2. If any discrepancies are found, the Region's designated representative may contact the laboratory to obtain additional information that could resolve any differences. If a discrepancy remains unresolved, the reviewer must use professional judgment to decide which value is the most accurate. Under these circumstances, the reviewer may determine that qualification of data is warranted. Note in the Data Review Narrative a description of the reasons for data qualification and the qualification that is applied to the data.
3. For non-aqueous samples, if the solids is less than 10.0%, use professional judgment for both detects and non-detects. If the percent solid for a soil sample is greater than or equal to 10.0% and less than 30.0%, use professional judgment to qualify detects and non-detects. If the percent solid for a soil sample is greater than or equal to 30.0%, detects and non-detects should not be qualified (see Table 11).
4. Note, for Contract Laboratory COR action, numerous or significant failures to accurately quantify the target compounds or to properly evaluate and adjust CRQLs.
5. Results between MDL and CRQL should be qualified as estimated "J".
6. Results < MDL should be reported at the CRQL and qualified "U". MDLs themselves should not be reported.

Table 11. Percent Solids Actions for Semivolatile Analysis for Non-Aqueous Samples

Criteria	Action	
	Detects	Non-detects
%Solids < 10.0%	Use professional judgment	Use professional judgment
10.0% ≤ %Solids ≤ 30.0%	Use professional judgment	Use professional judgment
%Solids > 30.0%	No qualification	No qualification

SAMPLE QUANTITATION

The sample quantitation evaluation is to verify laboratory quantitation results. In the space below, please show a minimum of one sample calculation:

Sample ID: JC19644-1 Analyte: 1,4-Dioxane RF: 0.619

$$\begin{aligned}
 [] &= (242629)(40)/(652029)(0.619) \\
 &= 24.04 \text{ ppm} \quad \text{Ok}
 \end{aligned}$$

DATA REVIEW WORKSHEETS

All criteria were met N/A
 Criteria were not met
 and/or see below

FIELD DUPLICATE PRECISION

Sample IDs: -

Matrix: -

Field duplicate samples may be taken and analyzed as an indication of overall precision. These analyses measure both field and lab precision; therefore, the results may have more variability than laboratory duplicates which only laboratory performance. It is also expected that soil duplicate results will have a greater variance than water matrices due to difficulties associated with collecting identical field duplicate samples.

The project QAPP should be reviewed for project-specific information.

Suggested criteria: if large RPD (> 50 %) is observed, confirm identification of the samples and note differences. If both samples and duplicate are <5 SQL, the RPD criteria is doubled.

COMPOUND	SQL ug/L	SAMPLE CONC.	DUPLICATE CONC.	RPD	ACTION
No field/laboratory duplicate analyzed as part of this data package. MS/MSD % and blank spike/blank spike duplicate recoveries RPD used to assess precision; RPD within the required criteria < 50 % for detected target analytes.					

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

OTHER ISSUES

A. System Performance

List samples qualified based on the degradation of system performance during sample analysis:

Sample ID	Comments	Actions
=====	=====	=====
_____	_____	_____
_____	_____	_____
_____	_____	_____
_____	_____	_____

Action:

Use professional judgment to qualify the data if it is determined that system performance has degraded during sample analyses. Inform the Contract Laboratory Program COR any action as a result of degradation of system performance which significantly affected the data.

B. Overall Assessment of Data

List samples qualified based on other issues:

Sample ID	Comments	Actions
=====	=====	=====
<u>No other issues that required the need to qualify the data. Results are valid and can be used for decision purposes.</u>		
_____	_____	_____

Action:

1. Use professional judgment to determine if there is any need to qualify data which were not qualified based on the Quality Control (QC) criteria previously discussed.
2. Write a brief narrative to give the user an indication of the analytical limitations of the data. Inform the Contract Laboratory COR the action, any inconsistency of the data with the Sample Delivery Group (SDG) Narrative. If sufficient information on the intended use and required quality of the data is available, the reviewer should include their assessment of the usability of the data within the given context. This may be used as part of a formal Data Quality Assessment (DQA).
3. Sometimes, due to dilutions, re-analysis or SIM/Scan runs are being performed, there will be multiple results for a single analyte from a single sample. The following criteria and professional judgment are used to determine which result should be reported:
 - The analysis with the lower CRQL
 - The analysis with the better QC results
 - The analysis with the higher results

EXECUTIVE NARRATIVE

SDG No: **JC19644** Laboratory: **Accutest, New Jersey**
Analysis: **SW846-8081B** Number of Samples: **5**

Location: **BMSMC, Building 5 Area**
Humacao, PR

SUMMARY: Five (5) samples were analyzed for selected pesticides following method SW846-8081B. The sample results were assessed according to USEPA data validation guidance documents in the following order of precedence *Hazardous Waste Support Section SOP No. HW-36A, Revision 0, June, 2015. SOM02.2. Pesticide Data Validation*. The QC criteria and data validation actions listed on the data review worksheets are from the primary guidance document, unless otherwise noted.

Results are valid and can be used for decision making purposes.


Critical issues: **None**
Major issues: **None**
Minor issues: **None**

Critical findings: **None**
Major findings: **None**
Minor findings: **None**

COMMENTS: Results are valid and can be used for decision making purposes.

Reviewers Name: Rafael Infante
Chemist License 1888

Signature:

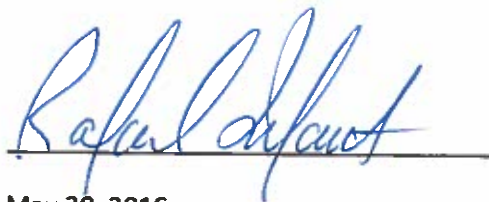A handwritten signature in blue ink, appearing to read 'Rafael Infante', is written over a horizontal line.

Date: May 30, 2016

SAMPLE ORGANIC DATA SAMPLE SUMMARY

Sample ID: JC19644-1
Sample location: BMSMC Building 5 Area
Sampling date: 3-May-16
Matrix: Groundwater

METHOD: 8081B

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Aldrin	0.0067	ug/L	1	-	U	Yes
alpha-BHC	0.0067	ug/L	1	-	U	Yes
beta-BHC	0.0067	ug/L	1	-	U	Yes
delta-BHC	0.0067	ug/L	1	-	U	Yes
gamma-BHC (Lindane)	0.0067	ug/L	1	-	U	Yes
alpha-Chlordane	0.0067	ug/L	1	-	U	Yes
gamma-Chlordane	0.0067	ug/L	1	-	U	Yes
Dieldrin	0.0067	ug/L	1	-	U	Yes
4,4'-DDD	0.0067	ug/L	1	-	U	Yes
4,4'-DDE	0.0067	ug/L	1	-	U	Yes
4,4'-DDT	0.0067	ug/L	1	-	U	Yes
Endrin	0.0067	ug/L	1	-	U	Yes
Endosulfan sulfate	0.0067	ug/L	1	-	U	Yes
Endrin aldehyde	0.0067	ug/L	1	-	U	Yes
Endrin ketone	0.0067	ug/L	1	-	U	Yes
Endosulfan-I	0.0067	ug/L	1	-	U	Yes
Endosulfan-II	0.0067	ug/L	1	-	U	Yes
Heptachlor	0.0067	ug/L	1	-	U	Yes
Heptachlor epoxide	0.0067	ug/L	1	-	U	Yes
Methoxychlor	0.013	ug/L	1	-	U	Yes
Toxaphene	0.17	ug/L	1	-	U	Yes

Sample ID: JC19644-2
Sample location: BMSMC Building 5 Area
Sampling date: 3-May-16
Matrix: Groundwater

METHOD: 8081B

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Aldrin	0.0067	ug/l	1	-	U	Yes
alpha-BHC	0.0067	ug/l	1	-	U	Yes
beta-BHC	0.0067	ug/l	1	-	U	Yes
delta-BHC	0.0067	ug/l	1	-	U	Yes
gamma-BHC (Lindane)	0.0067	ug/l	1	-	U	Yes
alpha-Chlordane	0.0067	ug/l	1	-	U	Yes
gamma-Chlordane	0.0067	ug/l	1	-	U	Yes
Dieldrin	0.0067	ug/l	1	-	U	Yes
4,4'-DDD	0.0067	ug/l	1	-	U	Yes
4,4'-DDE	0.0067	ug/l	1	-	U	Yes
4,4'-DDT	0.0067	ug/l	1	-	U	Yes
Endrin	0.0067	ug/l	1	-	U	Yes
Endosulfan sulfate	0.0067	ug/l	1	-	U	Yes
Endrin aldehyde	0.0067	ug/l	1	-	U	Yes
Endrin ketone	0.0067	ug/l	1	-	U	Yes
Endosulfan-I	0.0067	ug/l	1	-	U	Yes
Endosulfan-II	0.0067	ug/l	1	-	U	Yes
Heptachlor	0.0067	ug/l	1	-	U	Yes
Heptachlor epoxide	0.0067	ug/l	1	-	U	Yes
Methoxychlor	0.013	ug/l	1	-	U	Yes
Toxaphene	0.17	ug/l	1	-	U	Yes

Sample ID: JC19644-3
Sample location: BMSMC Building 5 Area
Sampling date: 3-May-16
Matrix: Groundwater

METHOD: 8081B

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Aldrin	0.0067	ug/l	1	-	U	Yes
alpha-BHC	0.0067	ug/l	1	-	U	Yes
beta-BHC	0.0067	ug/l	1	-	U	Yes
delta-BHC	0.0067	ug/l	1	-	U	Yes
gamma-BHC (Lindane)	0.0067	ug/l	1	-	U	Yes
alpha-Chlordane	0.0067	ug/l	1	-	U	Yes
gamma-Chlordane	0.0067	ug/l	1	-	U	Yes
Dieldrin	0.0067	ug/l	1	-	U	Yes
4,4'-DDD	0.0067	ug/l	1	-	U	Yes
4,4'-DDE	0.0067	ug/l	1	-	U	Yes
4,4'-DDT	0.0067	ug/l	1	-	U	Yes
Endrin	0.0067	ug/l	1	-	U	Yes
Endosulfan sulfate	0.0067	ug/l	1	-	U	Yes
Endrin aldehyde	0.0067	ug/l	1	-	U	Yes
Endrin ketone	0.0067	ug/l	1	-	U	Yes
Endosulfan-I	0.0067	ug/l	1	-	U	Yes
Endosulfan-II	0.0067	ug/l	1	-	U	Yes
Heptachlor	0.0067	ug/l	1	-	U	Yes
Heptachlor epoxide	0.0067	ug/l	1	-	U	Yes
Methoxychlor	0.013	ug/l	1	-	U	Yes
Toxaphene	0.17	ug/l	1	-	U	Yes

Sample ID: JC19644-4
Sample location: BMSMC Building 5 Area
Sampling date: 4-May-16
Matrix: Soil

METHOD: 8081B

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Aldrin	0.79	ug/kg		1 -	U	Yes
alpha-BHC	0.79	ug/kg	1	-	U	Yes
beta-BHC	0.79	ug/kg	1	-	U	Yes
delta-BHC	0.79	ug/kg	1	-	U	Yes
gamma-BHC (Lindane)	0.79	ug/kg	1	-	U	Yes
alpha-Chlordane	0.79	ug/kg	1	-	U	Yes
gamma-Chlordane	0.79	ug/kg	1	-	U	Yes
Dieldrin	0.79	ug/kg	1	-	U	Yes
4,4'-DDD	0.79	ug/kg	1	-	U	Yes
4,4'-DDE	0.79	ug/kg	1	-	U	Yes
4,4'-DDT	0.79	ug/kg	1	-	U	Yes
Endrin	0.79	ug/kg	1	-	U	Yes
Endosulfan sulfate	0.79	ug/kg	1	-	U	Yes
Endrin aldehyde	0.79	ug/kg	1	-	U	Yes
Endosulfan-I	0.79	ug/kg	1	-	U	Yes
Endosulfan-II	0.79	ug/kg	1	-	U	Yes
Heptachlor	0.79	ug/kg	1	-	U	Yes
Heptachlor epoxide	0.79	ug/kg	1	-	U	Yes
Methoxychlor	1.6	ug/kg	1	-	U	Yes
Endrin ketone	0.79	ug/kg	1	-	U	Yes
Toxaphene	20	ug/kg	1	-	U	Yes

Sample ID: JC19644-5
Sample location: BMSMC Building 5 Area
Sampling date: 4-May-16
Matrix: Soil

METHOD: 8081B

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Aldrin	0.83	ug/kg		1 -	U	Yes
alpha-BHC	0.83	ug/kg	1	-	U	Yes
beta-BHC	0.83	ug/kg	1	-	U	Yes
delta-BHC	0.83	ug/kg	1	-	U	Yes
gamma-BHC (Lindane)	0.83	ug/kg	1	-	U	Yes
alpha-Chlordane	0.83	ug/kg	1	-	U	Yes
gamma-Chlordane	0.83	ug/kg	1	-	U	Yes
Dieldrin	0.83	ug/kg	1	-	U	Yes
4,4'-DDD	0.83	ug/kg	1	-	U	Yes
4,4'-DDE	0.83	ug/kg	1	-	U	Yes
4,4'-DDT	0.83	ug/kg	1	-	U	Yes
Endrin	0.83	ug/kg	1	-	U	Yes
Endosulfan sulfate	0.83	ug/kg	1	-	U	Yes
Endrin aldehyde	0.83	ug/kg	1	-	U	Yes
Endosulfan-I	0.83	ug/kg	1	-	U	Yes
Endosulfan-II	0.83	ug/kg	1	-	U	Yes
Heptachlor	0.83	ug/kg	1	-	U	Yes
Heptachlor epoxide	0.83	ug/kg	1	-	U	Yes
Methoxychlor	1.7	ug/kg	1	-	U	Yes
Endrin ketone	0.83	ug/kg	1	-	U	Yes
Toxaphene	21	ug/kg	1	-	U	Yes

DATA REVIEW WORKSHEETS

Project/Case Number: JC19644Sampling Date: May_03-04, 2016Shipping Date: May_04, 2016EPA Region No.: 2

REVIEW OF PESTICIDE ORGANIC PACKAGE

The following guidelines for evaluating volatile organics were created to delineate required validation actions. This document will assist the reviewer in using professional judgment to make more informed decision and in better serving the needs of the data users. The sample results were assessed according to USEPA data validation guidance documents in the following order of precedence *Hazardous Waste Support Section SOP No. HW-36A, Revision 0, June, 2015. SOM02.2. Pesticide Data Validation*. The QC criteria and data validation actions listed on the data review worksheets are from the primary guidance document, unless otherwise noted.

The hardcopied (laboratory name) Accutest data package received has been reviewed and the quality control and performance data summarized. The data review for VOCs included:

Lab. Project/SDG No.: JC19644 Sample matrix: Groundwater
No. of Samples: 5

Trip blank No.: -
Field blank No.: -
Equipment blank No.: -
Field duplicate No.: -
Field spikes No.: -
QC audit samples: -

<input checked="" type="checkbox"/> Data Completeness	<input checked="" type="checkbox"/> Laboratory Control Spikes
<input checked="" type="checkbox"/> Holding Times	<input checked="" type="checkbox"/> Field Duplicates
<input type="checkbox"/> N/A GC/MS Tuning	<input checked="" type="checkbox"/> Calibrations
<input checked="" type="checkbox"/> Internal Standard Performance	<input checked="" type="checkbox"/> Compound Identifications
<input checked="" type="checkbox"/> Blanks	<input checked="" type="checkbox"/> Compound Quantitation
<input checked="" type="checkbox"/> Surrogate Recoveries	<input checked="" type="checkbox"/> Quantitation Limits
<input checked="" type="checkbox"/> Matrix Spike/Matrix Spike Duplicate	

Overall Comments: TCL_pesticides_list_by_SW846-8081B

Definition of Qualifiers:

J- Estimated results
U- Compound not detected
R- Rejected data
UJ- Estimated nondetect

Reviewer: Rafael Lopez
Date: May_30, 2016

DATA REVIEW WORKSHEETS

All criteria were met X
Criteria were not met
and/or see below

HOLDING TIMES

The objective of this parameter is to ascertain the validity of the results based on the holding time of the sample from time of collection to the time of analysis.

Complete table for all samples and note the analysis and/or preservation not within criteria

SAMPLE ID	DATE SAMPLED	DATE EXTRACTED/ANALYZED	ACTION
Samples properly preserved.			

Preservatives: All samples extracted and analyzed within the required criteria.

Criteria

Aqueous samples - seven (7) days from sample collection for extraction; 40 days from sample collection for analysis.

Non-aqueous samples – fourteen (14) days from sample collection for extraction; 40 days from sample collection for analysis.

Cooler temperature (Criteria: 4 ± 2 °C): 4°C - OK

Actions

Qualify aqueous sample results using preservation and technical holding time information as follows:

- If there is no evidence that the samples were properly preserved ($T = 4^{\circ}\text{C} \pm 2^{\circ}\text{C}$), and the samples were extracted or analyzed within the technical holding times, qualify detects as estimated (J) and non-detects as estimated (UJ).
- If there is no evidence that the samples were properly preserved ($T = 4^{\circ}\text{C} \pm 2^{\circ}\text{C}$), and the samples were extracted or analyzed outside the technical holding times, qualify detects as estimated (J) and non-detects as estimated (UJ).
- If the samples were properly preserved, and were extracted and analyzed within the technical holding times, no qualification of the data is necessary.
- If the samples were properly preserved, and were extracted or analyzed outside the technical holding times, qualify detects as estimated (J) and non-detects as estimated (UJ). Note in the Data Review Narrative that holding times were exceeded and the effect of exceeding the holding time on the resulting data.

DATA REVIEW WORKSHEETS

- e. Use professional judgment to qualify samples whose temperature upon receipt at the laboratory is either below 2 degrees centigrade or above 6 degrees centigrade.
- f. If technical holding times are grossly exceeded, use professional judgment to qualify the data.

Qualify non-aqueous sample results using preservation and technical holding time information as follows:

- a. If there is no evidence that the samples were properly preserved ($T = 4^{\circ}\text{C} \pm 2^{\circ}\text{C}$), and the samples were extracted or analyzed within the technical holding time, qualify detects as estimated (J) and non-detects as estimated (UJ).
- b. If there is no evidence that the samples were properly preserved ($T = 4^{\circ}\text{C} \pm 2^{\circ}\text{C}$), and the samples were extracted or analyzed outside the technical holding time, qualify detects as estimated (J) and non-detects as estimated (UJ).
- c. If the samples were properly preserved, and were extracted and analyzed within the technical holding time, no qualification of the data is necessary.
- d. If the samples were properly preserved, and were extracted or analyzed outside the technical holding time, qualify detects as estimated (J) and non-detects as estimated (UJ). Note in the Data Review Narrative that holding times were exceeded and the effect of exceeding the holding time on the resulting data.
- e. Use professional judgment to qualify samples whose temperature upon receipt at the laboratory is either below 2 degrees centigrade or above 6 degrees centigrade.
- f. If technical holding times are grossly exceeded, use professional judgment to qualify the data.

DATA REVIEW WORKSHEETS

All criteria were met X
Criteria were not met see below

GAS CHROMATOGRAPH WITH ELECTRON CAPTURE DETECTOR (GC/ECD) INSTRUMENT PERFORMANCE CHECK (SECTIONS 1 TO 5)

1. Resolution Check Mixture

Criteria

Is the resolution between two adjacent peaks in the Resolution Check Mixture C greater than or equal to 80.0% for all analytes for the primary column and greater than or equal to 50.0% for the confirmation column? Yes? or No?

Is the resolution between two adjacent peaks in the Resolution Check Mixture (A and B) greater than or equal to 60.0%? Yes? or No?

Note: If resolution criteria are not met, the quantitative results may not be accurate due to inadequate resolution. Qualitative identifications may also be questionable if coelution exists.

Action

- a. Qualify detects for target compounds that were not adequately resolved as tentatively identified (NJ).
- b. Qualify non-detected compounds as unusable (R).

2. Performance Evaluation Mixture (PEM) Resolution Criteria

Criteria

Is PEM analysis performed at the required frequency (at the end of each pesticide initial calibration sequence and every 12 hours)? Yes? or No?

Action

- a. If PEM is not performed at the required frequency, qualify all associated sample and blank results as unusable (R).

Criteria

Is PEM % Resolution < 90%? Yes? or No?

Action

- a. a. Qualify detects for target compounds that were not adequately resolved as tentatively identified (NJ).
- b. Qualify non-detected compounds as unusable (R).

DATA REVIEW WORKSHEETS

All criteria were met X
Criteria were not met see below

3. PEM 4,4'-DDT Breakdown

Criteria

Is the PEM 4,4'-DDT % Breakdown >20.0% and 4,4'-DDT is detected? Yes? or No?

Action

a. Qualify detects for 4,4'-DDT; detects for 4,4'-DDD; and detects for 4,4'-DDE as estimated (J)

Criteria

Is the PEM 4,4'-DDT % Breakdown >20.0% and 4,4'-DDT is not detected Yes? or No?

Action

- a. Qualify non-detects for 4,4'- DDT as unusable (R)
- b. Qualify detects for 4,4'-DDD as tentatively identified (NJ)
- c. Qualify detects for 4,4'-DDE as tentatively identified (NJ)

4. PEM Endrin Breakdown

Criteria

Is the PEM Endrin % Breakdown >20.0% and Endrin is detected? Yes? or No?

Action

a. Qualify detects for Endrin; detects for Endrin aldehyde; and detects for Endrin ketone as estimated (J)

Criteria

Is the PEM Endrin % Breakdown >20.0% and Endrin is not detected Yes? or No?

Action

- a. Qualify non-detects for Endrin as unusable (R)
- b. Qualify detects for Endrin aldehyde as tentatively identified (NJ)
- c. Qualify detects for Endrin ketone as tentatively identified (NJ)

All criteria were met X
Criteria were not met see below

5. Mid-point Individual Standard Mixture Resolution -

Criteria

Is the resolution between two adjacent peaks in the Resolution Check Mixture C greater than or equal to 80.0% for all analytes for the primary column and greater than or equal to 50.0% for the confirmation column?
Yes? or No?

Is the resolution between two adjacent peaks in the Resolution Check Mixture (A and B) greater than or equal to 90.0%?
Yes? or No?

Note: If resolution criteria are not met, the quantitative results may not be accurate due to inadequate resolution. Qualitative identifications may also be questionable if coelution exists.

Action

- a. Qualify detects for target compounds that were not adequately resolved as tentatively identified (NJ).
- b. Qualify non-detected compounds as unusable (R).

Criteria

Is mid-point individual standard mixture analysis performed at the required frequency (every 12 hours)?
Yes? or No?

Action

- a. If the mid-point individual standard mixture analysis is not performed at the required frequency, qualify all associated sample and blank results as unusable (R).

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below _____

CALIBRATION VERIFICATION

Compliance requirements for satisfactory instrument calibration are established to ensure that the instrument is capable of producing and maintaining acceptable quantitative data.

Date of initial calibration: 05/02/16 04/14/16
 Dates of initial calibration verification: 05/02/16 04/14/16
 Dates of continuing calibration: 05/06/16; 05/07/16 05/10/16
 Dates of final calibration: 05/06/16; 05/07/16 05/10/16
 Instrument ID numbers: GC4G HP_G1530_A
 Matrix/Level: Aqueous/low Aqueous/low

DATE	LAB ID#	FILE	CRITERIA OUT RFs, %RSD, %D, r	COMPOUND	SAMPLES AFFECTED
Initial and initial calibration verification within the guidance document performance criteria. Continuing calibration % differences meet the performance criteria in at least one of the column. Final calibration verification analyzed and results included in data package.					

Criteria

Are a five point calibration curve delivered with concentration levels as shown in Table 3 of SOP HW-36A, Revision 0, June, 2015? Yes? or No?

Actions

If the standard concentrations listed in Table 3 are not used, use professional judgment to evaluate the effect on the data

Criteria

Are RT Windows calculated correctly? Yes? or No?

Action

Recalculate the windows and use the corrected values for all evaluations.

Criteria

Are the Percent Relative Standard Deviation (%RSD) of the CFs for each of the single component target compounds less than or equal to 20.0%, except for alpha-BHC and delta-BHC?

Yes? or No?

DATA REVIEW WORKSHEETS

Are the %RSD of the CFs for alpha-BHC and delta-BHC less than or equal to 25.0%. Yes? or No?

Is the %RSD of the CFs for each of the Toxaphene peaks must be < 30% when 5-point ICAL is performed? Yes? or No?

Is the %RSD of the CFs for the two surrogates (tetrachloro-m-xylene and decachlorobiphenyl) less than or equal to 30.0%. Yes? or No?

Action

- If the %RSD criteria are not met, qualify detects as estimated (J) and use professional judgment to qualify non-detected target compounds.
- If the %RSD criteria are within allowable limits, no qualification of the data is necessary

Continuing Calibration Checks

Criteria

Is the continuing calibration standard analyzed at the acceptable time intervals? Yes? or No?

Action

- If more than 14 hours has elapsed from the injection of the instrument blank that begins an analytical sequence (opening CCV) and the injection of either a PEM or mid-point concentration of the Individual Standard Mixtures (A and B) or (C), qualify all data as unusable (R).
- If more than 12 hours has elapsed from the injection of the instrument blank that begins an analytical sequence (opening CCV) and the injection of the last sample or blank that is part of the same analytical sequence, qualify all data as unusable (R).
- If more than 72 hours has elapsed from the injection of the sample with a Toxaphene detection and the Toxaphene Calibration Verification Standard (CS3), qualify all data as unusable (R).

Criteria

Is the Percent Difference (%D) within $\pm 25.0\%$ for the PEM sample? Yes? or No?

Action

- Qualify associated detects as estimated (J) and non-detects as estimated (UJ).

Criteria

For the Calibration Verification Standard (CS3); is the Percent Difference (%D) within $\pm 25.0\%$? Yes? or No?

Action

Qualify associated detects as estimated (J) and non-detects as estimated (UJ).

DATA REVIEW WORKSHEETS

All criteria were met X
Criteria were not met
and/or see below

Criteria

Is the PEM 4,4'-DDT % Breakdown >20.0% and 4,4'-DDT is detected? Yes? or No?

Action

- a. Qualify detects for 4,4'-DDT; detects for 4,4'-DDD; and detects for 4,4'-DDE as estimated (J)
- b. Non-detected associated compounds are not qualified

Criteria

Is the PEM 4,4'-DDT % Breakdown >20.0% and 4,4'-DDT is not detected Yes? or No?

Action

- a. Qualify non-detects for 4,4'- DDT as unusable (R)
- b. Qualify detects for 4,4'-DDD as tentatively identified (NJ)
- c. Qualify detects for 4,4'-DDE as tentatively identified (NJ)

Criteria

Is the PEM Endrin % Breakdown >20.0% and Endrin is detected? Yes? or No?

Action

- a. Qualify detects for Endrin; detects for Endrin aldehyde; and detects for Endrin ketone as estimated (J)
- b. Non-detected associated compounds are not qualified

Criteria

Is the PEM Endrin % Breakdown >20.0% and Endrin is not detected Yes? or No?

Action

- a. Qualify non-detects for Endrin as unusable (R)
- b. Qualify detects for Endrin aldehyde as tentatively identified (NJ)
- c. Qualify detects for Endrin ketone as tentatively identified (NJ)

A separate worksheet should be filled for each initial curve

DATA REVIEW WORKSHEETS

All criteria were met X
Criteria were not met
and/or see below _____

BLANK ANALYSIS RESULTS (Sections 1 & 2)

The assessment of the blank analysis results is to determine the existence and magnitude of contamination problems. The criteria for evaluation of blanks apply only to blanks associated with the samples, including trip, equipment, and laboratory blanks. If problems with any blanks exist, all data associated with the case must be carefully evaluated to determine whether or not there is an inherent variability in the data for the case, or if the problem is an isolated occurrence not affecting other data.

List the contamination in the blanks below. High and low levels blanks must be treated separately.

CRQL concentration N/A

Laboratory blanks

DATE ANALYZED	LAB ID	LEVEL/ MATRIX	COMPOUND	CONCENTRATION UNITS
_No_target_analytes_detected_in_method_blanks_at_a_reporting_limit_of_0.01_and_0.001_ug/L.				

Field/Equipment/Trip blank

[illegible]

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

BLANK ANALYSIS RESULTS (Section 3)

Blank Actions

Action Levels (ALs) should be based upon the highest concentration of contaminant determined in any blank. Do not qualify any blank with another blank. The ALs for samples which have been diluted should be corrected for the sample dilution factor and/or % moisture, where applicable. No positive sample results should be reported unless the concentration of the compound in the samples exceeds the ALs:

The concentration of non-target compounds in all blanks must be less than or equal to 10 µg/L. The concentration of each target compound found in the method or field blanks must be less than its CRQL listed in the method.

Data concerning the field blanks are not evaluated as part of the CCS process. If field blanks are present, the data reviewer should evaluate this data in a similar fashion as the method blanks.

Specific actions are as follows:

Blank Actions for Pesticide Analyses

Blank Type	Blank Result	Sample Result	Action for Samples
Method, Sulfur Cleanup, Instrument, Field, TCLP/SPLP	Detects	Not detected	No qualification required
	< CRQL	< CRQL	Report CRQL value with a U
		≥ CRQL	No qualification required
	> CRQL	< CRQL	Report CRQL value with a U
		≥ CRQL and ≤ blank concentration	Report blank value for sample concentration with a U
		≥ CRQL and > blank concentration	No qualification required
	= CRQL	≤ CRQL	Report CRQL value with a U
		> CRQL	No qualification required
	Gross contamination	Detects	Report blank value for sample concentration with a U

DATA REVIEW WORKSHEETS

All criteria were met X
Criteria were not met
and/or see below _____

[illegible]

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

SURROGATE SPIKE RECOVERIES

Laboratory performance of individual samples is established by evaluation of surrogate spike recoveries. All samples are spiked with surrogate compounds prior to sample analysis. The accuracy of the analysis is measured by the surrogate percent recovery. Since the effects of the sample matrix are frequently outside the control of the laboratory and may present relatively unique problems, the validation of data is frequently subjective and demands analytical experience and professional judgment.

List the percent recoveries (%Rs) which do not meet the criteria for surrogate recovery.

Matrix: Aqueous

Lab Sample ID	Lab File ID	S1 a	S1 b	S2 a	S2 b
JC19644-1	4G67960.D	67	68	107	95
JC19644-2	4G67961.D	62	57	107	88
JC19644-3	4G67962.D	64	64	97	85
OP93665-BS1	4G67948.D	79	76	102	85
OP93665-MB1	4G67947.D	81	78	95	81
OP93665-MS	4G67996.D	65	66	96	90
OP93665-MSD	4G67997.D	65	69	94	90
Surrogate Compounds	Recovery Limits				

S1 = Tetrachloro-m-xylene 26-132% S2 = Decachlorobiphenyl 10-118%

(a) Recovery from GC signal #1

(b) Recovery from GC signal #2

Lab Sample ID	Lab File ID	S1 a	S1 b	S2 a	S2 b
JC19644-4	6G34906.D	109	107	106	120
JC19644-5	6G34907.D	104	107	109	119
OP93736-BS1	6G34898.D	109	108	107	122
OP93736-MB1	6G34897.D	89	89	88	100
OP93736-MS	6G34899.D	85	85	85	97
OP93736-MSD	6G34900.D	101	101	98	112
Surrogate Recovery Compounds	Recovery Limits				

S1 = Tetrachloro-m-xylene 24-136% S2 = Decachlorobiphenyl 10-153%

(a) Recovery from GC signal #1

(b) Recovery from GC signal #2

Note: Surrogate recoveries within laboratory control limits in project and QC samples.

DATA REVIEW WORKSHEETS

Actions:

- a. For any surrogate recovery greater than 150%, qualify detected target compounds as biased high (J+).
- b. Do not qualify non-detected target compounds for surrogate recovery > 150 %.
- c. If both surrogate recoveries are greater than or equal to 30% and less than or equal to 150%, no qualification of the data is necessary.
- d. For any surrogate recovery greater than or equal to 10% and less than 30%, qualify detected target compounds as biased low (J-).
- e. For any surrogate recovery greater than or equal to 10% and less than 30%, qualify non-detected target compounds as approximated (UJ).
- f. If low surrogate recoveries are from sample dilution, professional judgment should be used to determine if the resulting data should be qualified. If sample dilution is not a factor:
 - i. Qualify detected target compounds as biased low (J-).
 - ii. Qualify non-detected target compounds as unusable (R).
- g. If surrogate RTs in PEMs, Individual Standard Mixtures, samples, and blanks are outside of the RT Windows, the reviewer must use professional judgment to qualify data.
- h. If surrogate RTs are within RT windows, no qualification of the data is necessary.
- i. If the two surrogates were not added to all samples, MS/MSDs, standards, LCSs, and blanks, use professional judgment in qualifying data as missing surrogate analyte may not directly apply to target analytes.

Summary Surrogate Actions for Pesticide Analyses

Criteria	Action*	
	Detected Target Compounds	Non-detected Target Compounds
%R > 150%	J+	No qualification
30% < %R < 150%	No qualification	
10% < %R < 30%	J-	UJ
%R < 10% (sample dilution not a factor)	J-	R
%R < 10% (sample dilution is a factor)	Use professional judgment	
RT out of RT window	Use professional judgment	
RT within RT window	No qualification	

- * Use professional judgment in qualifying data, as surrogate recovery problems may not directly apply to target analytes.

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

MATRIX SPIKE/MATRIX SPIKE DUPLICATE (MS/MSD)

This data is generated to determine long term precision and accuracy in the analytical method for various matrices. This data alone cannot be used to evaluate the precision and accuracy of individual samples. If any % R in the MS or MSD falls outside the designated range, the reviewer should determine if there are matrix effects, i.e. LCS data are within the QC limits but MS/MSD data are outside QC limit.

1. MS/MSD Recoveries and Precision Criteria

Data for MS and MSDs will not be present unless requested by the Region.
 Notify the Contract Laboratory Program Project Officer (CLP PO) if a field blank was used for the MS and MSD, unless designated as such by the Region.

NOTE: For a Matrix Spike that does not meet criteria, apply the action to only the field sample used to prepare the Matrix Spike sample. If it is clearly stated in the data validation materials that the samples were taken through incremental sampling or some other method guaranteeing the homogeneity of the sample group, then the entire sample group may be qualified.

List the %Rs, RPD of the compounds which do not meet the criteria.

Sample ID: JC19644-1
 Sample ID: JC19823-1

Matrix/Level: Aqueous
 Matrix/Level: Soil

MS OR MSD	COMPOUND	% R	RPD	QC LIMITS	ACTION

Note: MS/MSD analyzed with this data package. % recoveries and RPD within laboratory control limits.

Action

No qualification of the data is necessary on MS and MSD data alone. However, using professional judgment, the validator may use the MS and MSD results in conjunction with other QC criteria and determine the need for some qualification of the data.

A separate worksheet should be used for each MS/MSD pair.

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below _____

LABORATORY CONTROL SAMPLE (LCS) ANALYSIS

This data is generated to determine accuracy of the analytical method for various matrices.

1. LCS Recoveries Criteria

LCS Spike Compound	Recovery Limits (%)
gamma-BHC	50 – 120
Heptachlor epoxide	50 – 150
Dieldrin	30 – 130
4,4'-DDE	50 – 150
Endrin	50 – 120
Endosulfan sulfate	50 – 120
trans-Chlordane	30 – 130
Tetrachloro-m-xylene (surrogate)	30 – 150
Decachlorobiphenyl (surrogate)	30 – 150

LCS concentrations: 0.167 ug/L; 16.7 ug/kg

List the %R of compounds which do not meet the criteria

LCS ID	COMPOUND	% R	QC LIMIT

Action

The following guidance is suggested for qualifying sample data for which the associated LCS does not meet the required criteria.

- If the LCS recovery exceeds the upper acceptance limit, qualify detected target compounds as estimated (J). Do not qualify non-detected target compounds.
- If the LCS recovery is less than the lower acceptance limit, qualify detected target compounds as estimated (J) and non-detects as unusable (R).
- Use professional judgment to qualify data for compounds other than those compounds that are included in the LCS.
- Use professional judgment to qualify non-LCS compounds. Take into account the compound class, compound recovery efficiency, analytical problems associated with each compound, and comparability in the performance of the LCS compound to the non-LCS compound.
- If the LCS recovery is within allowable limits, no qualification of the data is necessary.

DATA REVIEW WORKSHEETS

2. Frequency Criteria:

Where LCS analyzed at the required frequency and for each matrix? Yes or No.

If no, the data may be affected. Use professional judgment to determine the severity of the effect and qualify data accordingly. Discuss any actions below and list the samples affected.

Note: Blank spike/blank spike duplicate analyzed for aqueous matrices. % recoveries and RPD within laboratory control limits.

DATA REVIEW WORKSHEETS

All criteria were met _____
Criteria were not met _____
and/or see below N/A

FLORISIL CARTRIDGE PERFORMANCE CHECK

NOTE: Florisil cartridge cleanup is mandatory for all extracts.

Criteria

Is the Florisil cartridge performance check conducted at least once on each lot of cartridges used for sample cleanup or every 6 months, whichever is most frequent? Yes? or No? **N/A**

Criteria

Are the results for the Florisil Cartridge Performance Check solution included with the data package? Yes? or No? **N/A**

Note: If % criteria are not met, examine the raw data for the presence of polar interferences and use professional judgment in qualifying the data as follows:

Action:

- a. If the Percent Recovery is greater than 120% for any of the pesticide target compounds in the Florisil Cartridge Performance Check, qualify detected compounds as estimated (J). Do not qualify non-detected target compounds.
- b. If the Percent Recovery is greater than or equal to 80% and less than or equal to 120% for all the pesticide target compounds, no qualification of the data is necessary.
- c. If the Percent Recovery is greater than or equal to 10% and less than 80% for any of the pesticide target compounds in the Florisil Cartridge Performance Check, qualify detected target compounds as estimated (J) and non-detected target compounds as approximated (UJ).
- d. If the Percent Recovery is less than 10% for any of the pesticide target compounds in the Florisil Cartridge Performance Check, qualify detected compounds as estimated (J) and qualify non-detected target compounds as unusable (R).
- e. If the Percent Recovery of 2,4,5-trichlorophenol in the Florisil Cartridge Performance Check is greater than or equal to 5%, use professional judgment to qualify detected and non-detected target compounds, considering interference on the sample chromatogram.

Note: State in the Data Review Narrative potential effects on the sample data resulting from the Florisil Cartridge Performance Check analysis not yielding acceptable results.

Note: No information for florisil cartridge performance check included in data package. There is evidence that florisil cartridges were used in sample clean-up of soil samples. No action taken, professional judgment.

DATA REVIEW WORKSHEETS

All criteria were met N/A
Criteria were not met
and/or see below

GEL PERMEATION CHROMATOGRAPHY (GPC) PERFORMANCE CHECK

NOTE: GPC cleanup is mandatory for all soil samples.

If GPC criteria are not met, examine the raw data for the presence of high molecular weight contaminants; examine subsequent sample data for unusual peaks; and use professional judgment in qualifying the data. Notify the Contract Laboratory Program Project Officer (CLP PO) if the laboratory chooses to analyze samples under unacceptable GPC criteria.

Action:

- a. If the Percent Recovery is less than 10% for the pesticide compounds and surrogates during the GPC calibration check, the non-detected target compounds may be suspect, qualify detected compounds as estimated (J).
- b. If the Percent Recovery is less than 10% for the pesticide compounds and surrogates during the GPC calibration check, qualify all non-detected target compounds as unusable (R).
- c. If the Percent Recovery is greater than or equal to 10% and is less than 80% for any of the pesticide target compounds in the GPC calibration, qualify detected target compounds as estimated (J) and non-detected target compounds as approximated (UJ).
- d. If the Percent Recovery is greater than or equal to 80% and less than or equal to 120% for all the pesticide target compounds, no qualification of the data is necessary.
- e. If high recoveries (i.e., greater than 120%) were obtained for the pesticides and surrogates during the GPC calibration check, qualify detected compounds as estimated (J). Do not qualify non-detected target compounds.

Note: State in the Data Review Narrative potential effects on the sample data resulting from the GPC cleanup analyses not yielding acceptable results.

Note: No information for performance of GPC cleanup included in data package. No qualification of the data performed, professional judgment.

DATA REVIEW WORKSHEETS

All criteria were met X
Criteria were not met
and/or see below _____

TARGET COMPOUND IDENTIFICATION

Criteria:

1. Is Retention Times (RTs) of both of the surrogates and reported target compounds in each sample within the calculated RT Windows on both columns? Yes? or No?

2. Is the Tetrachloro-m-xylene (TCX) RT ± 0.05 minutes of the Mean RT (RT) determined from the initial calibration and Decachlorobiphenyl (DCB) within ± 0.10 minutes of the RT determined from the initial calibration? Yes? or No?

3. Is the Percent Difference (%D) for the detected mean concentrations of a pesticide target compound between the two Gas Chromatograph (GC) columns within the inclusive range of ± 25.0 %? Yes? or No?

4. When no analytes are identified in a sample; are the chromatograms from the analyses of the sample extract and the low-point standard of the initial calibration associated with those analyses on the same scaling factor? Yes? or No?

5. Does the chromatograms display the Single Component Pesticides (SCPs) detected in the sample and the largest peak of any multi-component analyte detected in the sample at less than full scale. Yes? or No?

6. If an extract is diluted; does the chromatogram display SCPs peaks between 10-100% of full scale, and multi-component analytes between 25-100% of full scale? Yes? or No? N/A

7. For any sample; does the baseline of the chromatogram return to below 50% of full scale before the elution time of alpha-BHC, and also return to below 25% of full scale after the elution time of alpha-BHC and before the elution time of DCB? Yes? or No?

8. If a chromatogram is replotted electronically to meet these requirements; is the scaling factor used displayed on the chromatogram, and both the initial chromatogram and the replotted chromatogram submitted in the data package. Yes? or No?

Action:

a. If the qualitative criteria for both columns were not met, all target compounds that are reported as detected should be considered non-detected.

b. Use professional judgment to assign an appropriate quantitation limit using the following guidance:

- i. If the detected target compound peak was sufficiently outside the pesticide RT Window, the reported values may be a false positive and should be replaced with the sample Contract Required Quantitation Limits (CRQL) value.

DATA REVIEW WORKSHEETS

- ii. If the detected target compound peak poses an interference with potential detection of another target peak, the reported value should be considered and qualified as unusable (R).
- c. If the data reviewer identifies a peak in both GC column analyses that falls within the appropriate RT Windows, but was reported as a non-detect, the compound may be a false negative. Use professional judgment to decide if the compound should be included.

Note: State in the Data Review Narrative all conclusions made regarding target compound identification.

- d. If the Toxaphene peak RT windows determined from the calibration overlap with SCPs or chromatographic interferences, use professional judgment to qualify the data.
- e. If target compounds were detected on both GC columns, and the Percent Difference between the two results is greater than 25.0%, consider the potential for coelution and use professional judgment to decide whether a much larger concentration obtained on one column versus the other indicates the presence of an interfering compound. If an interfering compound is indicated, use professional judgment to determine how best to report, and if necessary, qualify the data according to these guidelines.
- f. If Toxaphene exhibits a marginal pattern-matching quality, use professional judgment to establish whether the differences are due to environmental "weathering" (i.e., degradation of the earlier eluting peaks relative to the later eluting peaks). If the presence of Toxaphene is strongly suggested, report results as presumptively present (N).

GAS CHROMATOGRAPH/MASS SPECTROMETER (GC/MS) CONFIRMATION

NOTE: This confirmation is not usually provided by the laboratory. In cases where it is provided, use professional judgment to determine if data qualified with "C" can be salvaged if it was previously qualified as unusable (R).

Action:

- a. If the quantitative criteria for both columns were met (≥ 5.0 ng/ μ L for SCPs and ≥ 125 ng/ μ L for Toxaphene), determine whether GC/MS confirmation was performed. If it was performed, qualify the data using the following guidance:
 - i. If GC/MS confirmation was not required because the quantitative criteria for both columns was not met, but it was still performed, use professional judgment when evaluating the data to decide whether the detect should be qualified with "C".
 - ii. If GC/MS confirmation was performed, but unsuccessful for a target compound detected by GC/ECD analysis, qualify those detects as "X".

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

COMPOUND QUANTITATION AND REPORTED CONTRACT REQUIRED QUANTITATION LIMITS (CRQLS)

The sample quantitation evaluation is to verify laboratory quantitation results. In the space below, please show a minimum of one sample calculation:

Blank Spike alpha-BHC RF = 1.180

$$\begin{aligned} [] &= (253.9 \times 10^6)(50)/(439.8 \times 10^6)(1.180) \\ &= 24.5 \text{ ppb} \quad \text{Ok} \end{aligned}$$

Action:

- If sample quantitation is different from the reported value, qualify result as unusable (R).
- When a sample is analyzed at more than one dilution, the lowest CRQLs are used unless a QC exceedance dictates the use of the higher CRQLs from the diluted sample.
- Replace concentrations that exceed the calibration range in the original analysis by crossing out the "E" and its corresponding value on the original reporting form and substituting the data from the diluted sample.
- Results between the MDL and CRQL should be qualified as estimated (J).
- Results less than the MDL should be reported at the CRQL and qualified (U). MDLs themselves are not reported.
- For non-aqueous samples, if the percent moisture is less than 70.0%, no qualification of the data is necessary. If the percent moisture is greater than or equal to 70.0% and less than 90.0%, qualify detects as estimated (J) and non-detects as approximated (UJ). If the percent moisture is greater than or equal to 90.0%, qualify detects as estimated (J) and non-detects as unusable (R) (see Table).

Percent Moisture Actions for Pesticide Analysis for Non-Aqueous Samples

Criteria	Action	
	Detected Associated Compounds	Non-detected Associated Compounds
% Moisture < 70.0	No qualification	
70.0 < % Moisture < 90.0	J	UJ
% Moisture > 90.0	J	R

List samples which have $\leq 50\%$ solids

List samples which have $\leq 50\%$ solids

Dilution performed

[illegible]

DATA REVIEW WORKSHEETS

All criteria were met N/A
 Criteria were not met
 and/or see below

FIELD DUPLICATE PRECISION

NOTE: In the absence of QAPP guidance for validating data from field duplicates, the following action will be taken.

Field duplicates samples may be taken and analyzed as an indication of overall precision. These analyses measure both field and lab precision; therefore, the results may have more variability than laboratory duplicates which only laboratory performance. It is also expected that soil duplicate results will have a greater variance than water matrices due to difficulties associated with collecting identical field duplicate samples. Identify which samples within the data package are field duplicates. Estimate the relative percent difference (RPD) between the values for each compound. If large RPDs (> 50%) is observed, confirm identification of samples and note difference in the executive summary.

Sample IDs: -

Matrix: -

COMPOUND	SQL ug/L	SAMPLE CONC.	DUPLICATE CONC.	RPD	ACTION
No field/laboratory duplicate analyzed with this data package. LCS/LCSD % recoveries RPD used to assess precision. RPD within the required criteria of < 50 %.					

Actions:

a. Qualify as estimated positive results (J) and nondetects (UJ) for the compound that exceeded the above criteria. For organics, only the sample and duplicate will be qualified.

b. If an RPD cannot be calculated because one or both of the sample results is not detected, the following actions apply:

- i. If one sample result is not detected and the other is greater than 5x the SQL qualify (J/UJ).
- ii. If one sample value is not detected and the other is greater than 5x the SQL and the SQLs for the sample and duplicate are significantly different, use professional judgment to determine if qualification is appropriate.
- iii. If one sample value is not detected and the other is less than 5x, use professional judgment to determine if qualification is appropriate.
- iv. If both sample and duplicate results are not detected, no action is needed.

DATA REVIEW WORKSHEETS

OVERALL ASSESSMENT OF DATA

Action:

1. Use professional judgment to determine if there is any need to qualify data which were not qualified based on the Quality Control (QC) criteria previously discussed.
2. Write a brief narrative to give the user an indication of the analytical limitations of the data.

Note: The Contract Laboratory Program Project Officer (CLP PO) must be informed if any inconsistency of the data with the Sample Delivery Group (SDG) Narrative. If sufficient information on the intended use and required quality of the data is available, the reviewer should include their assessment of the usability of the data within the given context. This may be used as part of a formal Data Quality Assessment (DQA).

Overall assessment of the data: Results are valid; the data can be used for decision making purposes.

EXECUTIVE NARRATIVE

SDG No: **JC19644** Laboratory: **Accutest, New Jersey**
Analysis: **SW846-8015C** Number of Samples: **5**
Location: **BMSMC, Building 5 Area**
Humacao, PR

SUMMARY: Five (5) samples were analyzed for the low molecular weight alcohols (LMWAs) list following method SW846-8015C. The sample results were assessed according to USEPA data validation guidance documents in the following order of precedence: "Test Methods for Evaluating Solid Waste, Physical/Chemical Methods SW-846 (Final Update III, December 1996)," specifically for Methods 8000/8015C are utilized. The QC criteria and data validation actions listed on the data review worksheets are from the primary guidance document, unless otherwise noted.

Results are valid and can be used for decision making purposes.

Critical issues: **None**
Major: **None**
Minor: **None**


Critical findings: **None**
Major findings: **None**
Minor findings:

1. Initial, continuing, and final calibration verifications meets method specific criteria in at least one of the column, second column used for confirmation only. n-Butanol % difference over the method and guidance document performance criteria in continuing calibration verification. Results are qualified in affected sample as estimated (J) or (UJ).
2. MS/MSD recovery outside laboratory control limits for isobutyl alcohol and n-butyl alcohol in sample JC19644-1MS/-1MSD. Positive results qualified as estimated (J), non-detects were not qualified.
3. Blank spike/blank spike duplicate for isopropyl alcohol and n-butyl alcohol outside the upper control limits. Results were reported from 1st Signal. %RSD of initial calibration on 2nd signal excess method criteria (20%) so using for confirmation only. High percent recoveries and no associated positive found in the QC batch, results were accepted.

COMMENTS: Results are valid and can be used for decision making purposes.

Reviewers Name: **Rafael Infante**
Chemist License 1888

Signature:

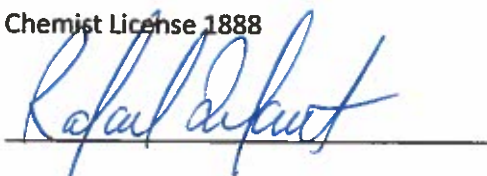

Date: **May 30, 2016**

SAMPLE ORGANIC DATA SAMPLE SUMMARY

Sample ID: JC19644-1
Sample location: BMSMC Building 5 Area
Sampling date: 5/3/2016
Matrix: Groundwater

METHOD: 8015C

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Ethanol	100	ug/l	1.0	-	U	Yes
Isobutyl Alcohol	100	ug/l	1.0	-	U	Yes
Isopropyl Alcohol	100	ug/l	1.0	-	U	Yes
n-Propyl Alcohol	100	ug/l	1.0	-	U	Yes
n-Butyl Alcohol	100	ug/l	1.0	-	UJ	Yes
sec-Butyl Alcohol	100	ug/l	1.0	-	U	Yes
Methanol	200	ug/l	1.0	-	U	Yes

Sample ID: JC19644-2
Sample location: BMSMC Building 5 Area
Sampling date: 5/3/2016
Matrix: Groundwater

METHOD: 8015C

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Ethanol	100	ug/l	1.0	-	U	Yes
Isobutyl Alcohol	100	ug/l	1.0	-	U	Yes
Isopropyl Alcohol	100	ug/l	1.0	-	U	Yes
n-Propyl Alcohol	100	ug/l	1.0	-	U	Yes
n-Butyl Alcohol	100	ug/l	1.0	-	UJ	Yes
sec-Butyl Alcohol	100	ug/l	1.0	-	U	Yes
Methanol	200	ug/l	1.0	-	U	Yes

Sample ID: JC19644-3
Sample location: BMSMC Building 5 Area
Sampling date: 5/3/2016
Matrix: Groundwater

METHOD: 8015C

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Ethanol	100	ug/l	1.0	-	U	Yes
Isobutyl Alcohol	100	ug/l	1.0	-	U	Yes
Isopropyl Alcohol	100	ug/l	1.0	-	U	Yes
n-Propyl Alcohol	100	ug/l	1.0	-	U	Yes
n-Butyl Alcohol	100	ug/l	1.0	-	UJ	Yes
sec-Butyl Alcohol	100	ug/l	1.0	-	U	Yes
Methanol	200	ug/l	1.0	-	U	Yes

Sample ID: JC19644-4
Sample location: BMSMC Building 5 Area
Sampling date: 5/4/2016
Matrix: Soil

METHOD: 8015C

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Ethanol	130	ug/kg	1.0	-	U	Yes
Isobutyl Alcohol	130	ug/kg	1.0	-	U	Yes
Isopropyl Alcohol	130	ug/kg	1.0	-	U	Yes
n-Propyl Alcohol	130	ug/kg	1.0	-	U	Yes
n-Butyl Alcohol	130	ug/kg	1.0	-	U	Yes
sec-Butyl Alcohol	130	ug/kg	1.0	-	U	Yes
Methanol	260	ug/kg	1.0	-	U	Yes

Sample ID: JC19644-5
Sample location: BMSMC Building 5 Area
Sampling date: 5/4/2016
Matrix: Soil

METHOD: 8015C

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Ethanol	130	ug/kg	1.0	-	U	Yes
Isobutyl Alcohol	130	ug/kg	1.0	-	U	Yes
Isopropyl Alcohol	130	ug/kg	1.0	-	U	Yes
n-Propyl Alcohol	130	ug/kg	1.0	-	U	Yes
n-Butyl Alcohol	130	ug/kg	1.0	-	U	Yes
sec-Butyl Alcohol	130	ug/kg	1.0	-	U	Yes
Methanol	250	ug/kg	1.0	-	U	Yes

DATA REVIEW WORKSHEETS

Project Number: JC19644
 Date: 05/03-04/2016
 Shipping Date: 05/04/2016
 EPA Region: 2

REVIEW OF VOLATILE ORGANIC PACKAGE

The following guidelines for evaluating volatile organics were created to delineate required validation actions. This document will assist the reviewer in using professional judgment to make more informed decision and in better serving the needs of the data users. The sample results were assessed according to USEPA data validation guidance documents in the following order of precedence: "Test Methods for Evaluating Solid Waste, Physical/Chemical Methods SW-846 (Final Update III, December 1996)," specifically for Methods 8000/8015C are utilized. The QC criteria and data validation actions listed on the data review worksheets are from the primary guidance document, unless otherwise noted.

The hardcopied (laboratory name) Accutest data package received has been reviewed and the quality control and performance data summarized. The modified data review for VOCs included:

Lab. Project/SDG No.: JC19644 Sample matrix: Soil/Groundwater
 No. of Samples: 5

Trip blank No.: -
 Field blank No.: -
 Equipment blank No.: -
 Field duplicate No.: -

<u>X</u> Data Completeness	<u>X</u> Laboratory Control Spikes
<u>X</u> Holding Times	<u>X</u> Field Duplicates
<u>N/A</u> GC/MS Tuning	<u>X</u> Calibrations
<u>N/A</u> Internal Standard Performance	<u>X</u> Compound Identifications
<u>X</u> Blanks	<u>X</u> Compound Quantitation
<u>X</u> Surrogate Recoveries	<u>X</u> Quantitation Limits
<u>X</u> Matrix Spike/Matrix Spike Duplicate	

Overall Comments: Low molecular weight alcohols by SW-846_8015C

Definition of Qualifiers:

J- Estimated results
 U- Compound not detected
 R- Rejected data
 UJ- Estimated nondetect

Reviewer: Rafael Infante
 Date: May 30, 2016

DATA REVIEW WORKSHEETS

DATA COMPLETENESS

MISSING INFORMATION

DATE LAB. CONTACTED

DATE RECEIVED

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

HOLDING TIMES

The objective of this parameter is to ascertain the validity of the results based on the holding time of the sample from time of collection to the time of analysis.

Complete table for all samples and note the analysis and/or preservation not within criteria

SAMPLE ID	DATE SAMPLED	DATE ANALYZED	pH	ACTION
All samples analyzed within the recommended method holding time. All samples properly preserved.				

Criteria

Aqueous samples – 14 days from sample collection for preserved samples ($\text{pH} \leq 2$, 4°C), no air bubbles.

Aqueous samples – 7 days from sample collection for unpreserved samples, 4°C , no air bubbles.

Soil samples- 7 days from sample collection.

Cooler temperature (Criteria: $4 \pm 2^{\circ}\text{C}$): 4°C

Actions

If the VOCs vial(s) have air bubbles, estimate positive results (J) and reject nondetects (R).

If the % solids of soil samples is 10-50%, estimate positive results (J) and nondetects (UJ)

If the % solid of soil samples is $< 10\%$, estimate positive results (J) and reject nondetects (R).

If holding times are exceeded but < 14 days beyond criteria, estimate positive results (J) and nondetects (UJ).

If holding times are exceeded but < 28 days beyond criteria, estimate positive results (J) and reject nondetects (R).

If holding times are grossly exceeded (> 28 days beyond criteria), reject all results (R).

If samples were not iced or if the ice were melted ($> 10^{\circ}\text{C}$), estimate positive results (J) and nondetects (UJ).

DATA REVIEW WORKSHEETS

All criteria were met __N/A__
Criteria were not met see below

GC/MS TUNING

The assessment of the tuning results is to determine if the sample instrumentation is within the standard tuning QC limits

__N/A__ The BFB performance results were reviewed and found to be within the specified criteria.

__N/A__ BFB tuning was performed for every 12 hours of sample analysis.

If no, use professional judgment to determine whether the associated data should be accepted, qualified or rejected.

List the samples affected: _____

If mass calibration is in error, all associated data are rejected.

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below _____

CALIBRATION VERIFICATION

Compliance requirements for satisfactory instrument calibration are established to ensure that the instrument is capable of producing and maintaining acceptable quantitative data.

Date of initial calibration: 05/09/16; 05/17/16
 Dates of initial calibration verification: 05/09/16; 05/17/16
 Dates of continuing calibration verification: 05/16/16; 05/18/16
 Dates of final calibration verification: 05/16/16; 05/18/16
 Instrument ID number: GCGH
 Matrix/Level: Aqueous/low

DATE	LAB FILE ID#	CRITERIA OUT RFs, %RSD, %D, r	COMPOUND	SAMPLES AFFECTED
05/16/16	CC5269-5000	-24.9/61.0	1-Butanol	JC19644-1; -2; -3; - 1MS-1MSD

Note: Initial, continuing, and final calibration verifications meets method specific criteria except in the cases described in this document. Results are qualified in affected sample as estimated (J) or (UJ).

Criteria

All RFs must be > 0.05 regardless of method requirements for SPCC.

All %RSD must be $\leq 15\%$ regardless of method requirements for CCC.

All %Ds must be $\leq 20\%$ regardless of method requirements for CCC.

It should be noted that Region 2 SOP HW-24 does not specify criterion for the curve correlation coefficient (r). A limit for r of ≥ 0.995 has therefore been utilized as professional judgment.

Actions

If any compound has an initial RF or a continuing RF of < 0.05 , estimate positive results (J) and reject nondetects (R), regardless of method requirements.

If any compound has a %RSD $> 15\%$, estimate positive results (J) and use professional judgment to qualify nondetects.

If any compound has a %RSD $> 90\%$, estimate positive results (J) and reject nondetects (R).

If any compound has a % D $> 20\%$, estimate positive results (J) and reject nondetects (R).

If any compound has a % D $> 20\%$, estimate positive results (J) and nondetects (UJ).

If any compound has a % D $> 90\%$, estimate positive results (J) and reject nondetects (R).

If any compound has $r < 0.995$, estimate positive results and nondetects.

A separate worksheet should be filled for each initial curve

5

5

5

5

5

5

5

5

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

V B. BLANK ANALYSIS RESULTS (Section 3)

Blank Actions

Action Levels (ALs) should be based upon the highest concentration of contaminant determined in any blank. Do not qualify any blank with another blank. The ALs for samples which have been diluted should be corrected for the sample dilution factor and/or % moisture, where applicable. No positive sample results should be reported unless the concentration of the compound in the samples exceeds the ALs:

ALs = 10x the amount of common contaminants (methylene chloride, acetone, 2-butanone, and toluene)
 ALs = 5x for any other compounds

Specific actions are as follows:

If the concentration is < sample quantitation limit (SQL) and \leq AL, report the compound as not detected (U) at the SQL.

If the concentration is \geq SQL but \leq AL, report the compound as not detected (U) at the reported concentration.

If the concentration is \geq SQL and > AL, report the concentration unqualified.

Notes:

High and low level blanks must be treated separately

Compounds qualified "U" for blank contamination are still considered "hits" when qualifying for calibration criteria.

CONTAMINATION SOURCE/LEVEL	COMPOUND	CONC/UNITS	AL/UNITS	SQL	AFFECTED SAMPLES

DATA REVIEW WORKSHEETS

All criteria were met ☒
 Criteria were not met
 and/or see below _____

SURROGATE SPIKE RECOVERIES

Laboratory performance of individual samples is established by evaluation of surrogate spike recoveries. All samples are spiked with surrogate compounds prior to sample analysis. The accuracy of the analysis is measured by the surrogate percent recovery. Since the effects of the sample matrix are frequently outside the control of the laboratory and may present relatively unique problems, the validation of data is frequently subjective and demands analytical experience and professional judgment.

List the percent recoveries (%Rs) which do not meet the criteria for surrogate recovery.

Matrix: solid/aqueous

SAMPLE ID	SURROGATE COMPOUND				ACTION
	Hexanol	DBFM	TOL-d8	BFB	

All surrogate recoveries within laboratory control limits. _____

QC Limits* (Aqueous)

_____ LL to UL _____ 73 to 123 _____ to _____ to _____ to _____

QC Limits* (Solid-Low)

_____ LL to UL _____ 69 to 121 _____ to _____ to _____ to _____

QC Limits* (Solid-Med)

_____ LL to UL _____ to _____ to _____ to _____ to _____

1,2-DCA = 1,2-Dichloromethane-d4

TOL-d8 = Toluene-d8

DBFM = Dibromofluoromethane

BFB = Bromofluorobenzene

* QC limits are laboratory in-house performance criteria, LL = lower limit, UL = upper limit.

* If QC limits are not available, use limits of 80 – 120 % for aqueous and 70 – 130 % for solid samples.

Actions:

QUALITY	%R < 10%	%R = 10% - LL	%R > UL
Positive results	J	J	J
Nondetects results	R	UJ	Accept

Surrogate action should be applied:

If one or more surrogate in the VOC fraction is out of specification, but has a recovery of > 10%.

If any one surrogate in a fraction shows < 10 % recovery.

DATA REVIEW WORKSHEETS

All criteria were met ☒
 Criteria were not met ☐
 and/or see below ☐

VII. A MATRIX SPIKE/MATRIX SPIKE DUPLICATE (MS/MSD)

This data is generated to determine long term precision and accuracy in the analytical method for various matrices. This data alone cannot be used to evaluate the precision and accuracy of individual samples. If any % R in the MS or MSD falls outside the designated range, the reviewer should determine if there are matrix effects, i.e. LCS data are within the QC limits but MS/MSD data are outside QC limit.

1. MS/MSD Recoveries and Precision Criteria

The laboratory should use one MS and a duplicate analysis of an unspiked field sample if target analytes are expected in the sample. If target analytes are not expected, MS/MSD should be analyzed.

List the %Rs, RPD of the compounds which do not meet the criteria.

Sample ID: JC19644-1MS/-MSD Matrix/Level: Aqueous
 Sample ID: JC19914-2MS/-MSD Matrix/Level: Soil

MS OR MSD	COMPOUND	% R	RPD	QC LIMITS	ACTION
-----------	----------	-----	-----	-----------	--------

MS/MSD%_recoveries_and_RPD_within_laboratory_control_limits_except_for_the_followings:_

The QC reported here applies to the following samples: Method: SW846 8015C MOD
 JC19644-1; JC19644-2; JC19644-3

Compound	JC19644-1 mg/l Q	Spike mg/l	MS mg/l	MS %	Spike mg/l	MSD mg/l	MSD %	RPD	Limits Rec/RPD
Isobutyl Alcohol	ND	5000	6360	127	5000	6620	132*	4	69-131/25
n-Butyl Alcohol	ND	5000	6710	134* a	5000	6890	138* a	3	63-131/25

Note: Positive results are qualified as estimated in sample JC19644-1, non-detects results are accepted.

- * QC limits are laboratory in-house performance criteria, LL = lower limit, UL = upper limit.
- * If QC limits are not available, use limits of 70 – 130 %.

Actions:

QUALITY	%R < LL	%R > UL
Positive results	J	J
Nondetects results	R	Accept

DATA REVIEW WORKSHEETS

All criteria were met X
 Criteria were not met
 and/or see below

MS/MSD criteria apply only to the unspiked sample, its dilutions, and the associated MS/MSD samples:

If the % R for the affected compounds were < LL (or 70 %), qualify positive results (J) and nondetects (UJ).

If the % R for the affected compounds were > UL (or 130 %), only qualify positive results (J).

If 25 % or more of all MS/MSD %R were < LL (or 70 %) or if two or more MS/MSD %Rs were < 10%, qualify all positive results (J) and reject nondetects (R).

VII. B MATRIX SPIKE/MATRIX SPIKE DUPLICATE

MS/MSD – Unspiked Compounds

It should be noted that Region 2 SOP HW-24 does not specify a MS/MSD criteria for the unspiked compounds in the sample. A %RSD of < 50% has therefore been utilized as professional judgment.

If all target analytes were spiked in the MS/MSD, this review element is not applicable.

List the %RSD of the compounds which do not meet the criteria.

Sample ID: _____ - _____ Matrix/Level/Unit: _____ - _____

COMPOUND	SAMPLE CONC.	MS CONC.	MSD CONC.	% RSD	ACTION

Actions:

* If the % RSD > 50, qualify the positive result in the unspiked samples as estimated (J).

* If the % RSD is not calculated (NC) due to nondetected value, use professional judgment to qualify the data.

A separate worksheet should be used for each MS/MSD pair.

DATA REVIEW WORKSHEETS

All criteria were met X
Criteria were not met
and/or see below _____

VIII. LABORATORY CONTROL SAMPLE (LCS) ANALYSIS

This data is generated to determine accuracy of the analytical method for various matrices.

1. LCS Recoveries Criteria

Where LCS spiked with the same analyte at the same concentrations as the MS/MSD? **Yes**
or No. If no make note in data review memo.

List the %R of compounds which do not meet the criteria

LCS ID	COMPOUND	% R	QC LIMIT
____ Recoveries within laboratory control limits except for the following: _____			
__GGH5285-BS a_(Aqueous)	Isopropyl alcohol	174*_b	76_-_121
	n-Butyl alcohol	126*_b	67_-_116

(a) Reported from 1st Signal. %RSD of initial calibration on 2nd signal excess method criteria (20%) so using for confirmation only.

(b) High percent recoveries and no associated positive found in the QC batch.

Note: No action taken, analytes not detected in samples.

- * QC limits are laboratory in-house performance criteria, LL = lower limit, UL = upper limit.
- * If QC limits are not available, use limits of 70 – 130 %.

Actions:

QUALITY	%R < LL	%R > UL
Positive results	J	J
Nondetects results	R	Accept

All analytes in the associated sample results are qualified for the following criteria.

If 25 % of the LCS recoveries were < LL (or 70 %), qualify all positive results (j) and reject nondetects (R).

If two or more LCS were below 10 %, qualify all positive results as (J) and reject nondetects (R).

2. Frequency Criteria:

Where LCS analyzed at the required frequency and for each matrix? Yes or No.

If no, the data may be affected. Use professional judgment to determine the severity of the effect and qualify data accordingly. Discuss any actions below and list the samples affected.

DATA REVIEW WORKSHEETS

All criteria were met N/A
 Criteria were not met
 and/or see below

IX. FIELD/LABORATORY DUPLICATE PRECISION

Sample IDs: - Matrix: -

Field/laboratory duplicates samples may be taken and analyzed as an indication of overall precision. These analyses measure both field and lab precision; therefore, the results may have more variability than laboratory duplicates which only laboratory performance. It is also expected that soil duplicate results will have a greater variance than water matrices due to difficulties associated with collecting identical field duplicate samples.

The project QAPP should be reviewed for project-specific information.

Suggested criteria: RPD \pm 30% for aqueous samples, RPD \pm 50 % for solid samples. If both samples and duplicate are <5 SQL, the RPD criteria is doubled.

COMPOUND	SQL	SAMPLE CONC.	DUPLICATE CONC.	RPD	ACTION
No field/laboratory duplicate analyzed with this data package. MS/MSD % recoveries RPD used to assess precision. RPD within laboratory and generally acceptable control limits.					

Actions:

Qualify as estimated positive results (J) and nondetects (UJ) for the compound that exceeded the above criteria. For organics, only the sample and duplicate will be qualified.

If an RPD cannot be calculated because one or both of the sample results is not detected, the following actions apply:

If one sample result is not detected and the other is greater than 5x the SQL qualify (J/UJ).

If one sample value is not detected and the other is greater than 5x the SQL and the SQLs for the sample and duplicate are significantly different, use professional judgment to determine if qualification is appropriate.

If one sample value is not detected and the other is less than 5x, use professional judgment to determine if qualification is appropriate.

If both sample and duplicate results are not detected, no action is needed.

All criteria were met N/A
Criteria were not met
and/or see below _____

X. INTERNAL STANDARD PERFORMANCE

The assessment of the internal standard (IS) parameter is used to assist the data reviewer in determining the condition of the analytical instrumentation.

List the internal standard area of samples which do not meet the criteria.

- * Area of +100% or -50% of the IS area in the associated calibration standard.
* Retention time (RT) within 30 seconds of the IS area in the associated calibration standard.

DATE	SAMPLE ID	IS OUT	IS AREA	ACCEPTABLE RANGE	ACTION
------	-----------	--------	---------	------------------	--------

[illegible]

Actions:

1. IS actions should be applied to the compound quantitated with the out-of-control ISs

QUALITY	IS AREA < -25%	IS AREA = -25 % TO - 50%	IS AREA > + 100%
Positive results	J	J	J
Nondetected results	R	UJ	ACCEPT

2. If a IS retention time varies more than 30 seconds, the chromatographic profile for that sample must be examined to determine if any false positive or negative exists. For shifts of a large magnitude, the reviewer may consider partial or total rejection of the data for the sample fraction.

DATA REVIEW WORKSHEETS

All criteria were met X
Criteria were not met
and/or see below

XII. SAMPLE QUANTITATION

The sample quantitation evaluation is to verify laboratory quantitation results. In the space below, please show a minimum of one sample calculation:

JC19644-1

Hexanol

RF = 76.51

$$[] = (398784)/(76.51)$$

$$= 5,212 \text{ ppm OK}$$

— — — — —

All criteria were met X
Criteria were not met
and/or see below _____

XII. QUANTITATION LIMITS

A. Dilution performed

[illegible]

B. Percent Solids

List samples which have $\leq 50\%$ solids

Actions:

If the % solids of a soil sample is 10-50%, estimate positive results (J) and nondetects (UJ)

If the % solids of a soil sample is $< 10\%$, estimate positive results (J) and reject nondetects (R)