CETIFICATION

SDG No:

JC21103

Laboratory:

Accutest, New Jersey

Site:

BMS, Building 5 Area, PR

Matrix:

Soil/Groundwater

Humacao, PR

SUMMARY:

Groundwater and soil samples (Table 1) were collected on the BMSMC facility – Building 5 Area. The BMSMC facility is located in Humacao, PR. Samples were taken May 25-26, 2016 and were analyzed in Accutest Laboratory of Dayton, New Jersey for the ABN TCL Special List (1,4-Dioxane and Naphthalene were analyzed following the SIM technique); TCL pesticides list; and for low molecular weight alcohols (LMWA) the results were reported under SDG No.: JC21103. Results were validated using the latest validation guidelines (July, 2015) of the EPA Hazardous Waste Support Section. The analyses performed are shown in Table 1. Individual data review worksheets are enclosed for each target analyte group. The data sample organic data samples summary form shows for analytes results that were qualified.

In summary the results are valid and can be used for decision taking purposes.

Table 1. Samples analyzed and analysis performed

SAMPLE ID	SAMPLE DESCRIPTION	MATRIX	ANALYSIS PERFORMED
JC21103-1	RA7-GWS	Groundwater	ABN TCL special list; 1,-4- dioxane and Naphthalene (SIM); Pesticides TCL list; LMWA
JC21103-2	EB05252016	Equipment Blank	ABN TCL special list; 1,-4- dioxane and Naphthalene (SIM); Pesticides TCL list; LMWA
JC21103-3	SB103 (2.5-3.5)	Soil	ABN TCL special list; 1,-4- dioxane and Naphthalene (SIM); LMWA
JC21103-4	SB103D (2.5-3.5)	Soil	ABN TCL special list; 1,-4- dioxane and Naphthalene (SIM); LMWA
JC21103-5	SB103 (6-7)	Soil	ABN TCL special list; 1,-4- dioxane and Naphthalene (SIM); LMWA
JC21103-6	MVV-235 (5-6)	Soil	ABN TCL special list; 1,-4- dioxane and Naphthalene (SIM); LMWA
JC21103-7	SB103-GWD	Groundwater	ABN TCL special list; 1,-4- dioxane and Naphthalene (SIM); Pesticides TCL list; LMWA
JC21103-8	SB103-GWS	Groundwater	ABN TCL special list; 1,-4- dioxane and Naphthalene (SIM); Pesticides TCL list; LMWA

Reviewer Name:

Rafael Infante

Chemist License 1888

Signature:

Date:

June 22, 2016

STOCK STOCK

Page 1 of 3

Client Sample ID: Lab Sample ID:

RA7-GWS JC21103-1

Matrix:

AQ - Ground Water

SW846 8270D SW846 3510C

Date Sampled: Date Received:

05/25/16 05/27/16

Method:

Percent Solids:

Project:

BMSMC, Building 5 Area, PR

Prep Batch

Run #1

Z111120.D

File ID

DF 1

Analyzed By 05/28/16 IJ

Prep Date 05/27/16

OP94327

Q

Analytical Batch EZ5550

Run #2

Final Volume Initial Volume

Run #1

920 ml

1.0 ml

Run #2

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units
95-57-8	2-Chlorophenol	ND	5.4	0.89	ug/l
59-50-7	4-Chloro-3-methyl phenol	ND	5.4	0.97	ug/l
120-83-2	2,4-Dichlorophenol	ND	2.2	1.4	ug/l
105-67-9	2,4-Dimethylphenol	ND	5.4	2.7	ug/l
51-28-5	2,4-Dinitrophenol	ND	11	1.7	ug/l
534-52-1	4,6-Dinitro-o-cresol	ND	5.4	1.4	ug/l
95-48-7	2-Methylphenol	ND	2.2	0.97	ug/l
	3&4-Methylphenol	ND	2.2	0.96	ug/l
88-75-5	2-Nitrophenol	ND	5.4	1.0	ug/l
100-02-7	4-Nitrophenol	ND	11	1.3	ug/l
87-86-5	Pentachlorophenol	ND	5.4	1.5	ug/l
108-95-2	Phenol	ND	2.2	0.43	ug/l
58-90-2	2,3,4,6-Tetrachlorophenol	ND	5.4	1.6	ug/l
95-95-4	2,4,5-Trichlorophenol	ND	5.4	1.4	ug/l
88-06-2	2,4,6-Trichlorophenol	ND.	5.4	1.0	ug/l
83-32-9	Acenaphthene	ND	1.1	0.21	ug/l
208-96-8	Acenaphthylene	ND	1.1	0.15	ug/l
98-86-2	Acetophenone	ND	2.2	0.23	ug/l
120-12-7	Anthracene	ND	1.1	0.23	ug/l
1912-24-9	Atrazine	ND	2.2	0.49	ug/l
100-52-7	Benzaldehyde	ND	5.4	0.31	ug/l
56-55-3	Benzo(a)anthracene	ND	1.1	0.22	ug/l
50-32-8	Benzo(a)pyrene	ND	1.1	0.23	ug/l
205-99-2	Benzo(b)fluoranthene	ND	1.1	0.22	ug/l
191-24-2	Benzo(g,h,i)perylene	ND	1.1	0.37	ug/l
207-08-9	Benzo(k)fluoranthene	ND	1.1	0.22	ug/l
101-55-3	4-Bromophenyl phenyl ether	ND	2.2	0.44	ug/l
85-68-7	Butyl benzyl phthalate	ND	2.2	0.50	ug/l
92-52-4	1,1'-Biphenyl	ND	1.1	0.23	ug/l
91-58-7	2-Chloronaphthalene	ND	2.2	0.26	ug/l
106-47-8	4-Chloroaniline	ND	5.4	0.37	ug/l
00.74.0	C11-	AID		0.05	

fuel Infanta Méndez

ND = Not detected

86-74-8

MDL = Method Detection Limit

ND

RL = Reporting Limit

E = Indicates value exceeds calibration range

Carbazole

J = Indicates an estimated value

ug/l

0.25

B = Indicates analyte found in associated method blank



Project:

Client Sample ID: **RA7-GWS** Lab Sample ID: JC21103-1

Matrix: AQ - Ground Water Method: SW846 8270D SW846 3510C

BMSMC, Building 5 Area, PR

Report of Analysis

Date Sampled: 05/25/16 Date Received: 05/27/16 Percent Solids: n/a

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units	Q	
105-60-2	Caprolactam	ND	2.2	0.71	ug/l		
218-01-9	Chrysene	ND	1.1	0.19	ug/l		
111-91-1	bis(2-Chloroethoxy)methane	ND	2.2	0.30	ug/l		
111-44-4	bis(2-Chloroethyl)ether	ND	2.2	0.27	ug/l		
108-60-1	bis(2-Chloroisopropyl)ether	ND	2.2	0.44	ug/l		
7005-72-3	4-Chlorophenyl phenyl ether	ND	2.2	0.40	ug/l		
121-14-2	2,4-Dinitrotoluene	ND	1.1	0.60	ug/I		
606-20-2	2,6-Dinitrotoluene	ND	1.1	0.52	ug/l		
91-94-1	3,3'-Dichlorobenzidine	ND	2.2	0.55	ug/l		
123-91-1	1,4-Dioxane	6.6	1.1	0.71	ug/l		
53-70-3	Dibenzo(a,h)anthracene	ND	1.1	0.36	ug/l		
132-64-9	Dibenzofuran	ND	5.4	0.24	ug/l		
84-74-2	Di-n-butyl phthalate	ND	2.2	0.54	ug/l		
117-84-0	Di-n-octyl phthalate	ND	2.2	0.25	ug/l		
84-66-2	Diethyl phthalate	ND	2.2	0.28	ug/l		
131-11-3	Dimethyl phthalate	ND	2.2	0.24	ug/l		
117-81-7	bis(2-Ethylhexyl)phthalate	ND	2.2	1.8	ug/l		
206-44-0	Fluoranthene	ND	1.1	0.18	ug/l		
86-73-7	Fluorene	ND	1.1	0.19	ug/l		
118-74-1	Hexachlorobenzene	ND	1.1	0.35	ug/l		
87-68-3	Hexachlorobutadiene	ND	1.1	0.53	ug/l		
77-47-4	Hexachlorocyclopentadiene	ND	11	3.0	ug/l		
67-72-1	Hexachloroethane	ND	2.2	0.42	ug/l		
193-39-5	Indeno(1,2,3-cd)pyrene	ND	1.1	0.36	ug/I		
78-59-1	Isophorone	ND	2.2	0.30	ug/l		
90-12-0	1-Methylnaphthalene	ND	1.1	0.29	ug/l		
91-57-6	2-Methylnaphthalene	ND	1.1	0.23	ug/l		
88-74-4	2-Nitroaniline	ND	5.4	0.30	ug/l		
99-09-2	3-Nitroaniline	ND	5.4	0.42	ug/l		
100-01-6	4-Nitroaniline	ND	5.4	0.48	ug/l		14
98-95-3	Nitrobenzene	ND	2.2	0.70	ug/l	/	gat
621-64-7	N-Nitroso-di-n-propylamine	ND	2.2	0.52	ug/l	/5	5%
86-30-6	N-Nitrosodiphenylamine	ND	5.4	0.24	ug/l	12	1/ 1
85-01-8	Phenanthrene	ND	1.1	0.19	ug/l	55.	15
129-00-0	Pyrene	ND	1.1	0.24	ug/l	/ •	1
95-94-3	1,2,4,5-Tetrachlorobenzene	ND	2.2	0.40	ug/l	EST, 17.3	11
CAS No.	Surrogate Recoveries	Run#1	Run# 2	Lim	its		177
367-12-4	2-Fluorophenol	42%		14-8	8%		



ND = Not detected MDL = Method Detection Limit

RL = Reporting Limit

E = Indicates value exceeds calibration range

J = Indicates an estimated value

B = Indicates analyte found in associated method blank

Client Sample ID: **RA7-GWS** Lab Sample ID: JC21103-1

Matrix: Method:

Project:

AQ - Ground Water

SW846 8270D SW846 3510C BMSMC, Building 5 Area, PR

Report of Analysis

Date Sampled: Date Received:

05/25/16 05/27/16

Percent Solids: n/a

ABN TCL Special List

CAS No.	Surrogate Recoveries	Run#1	Run# 2	Limits
4165-62-2	Phenol-d5	32%		10-110%
118-79-6	2,4,6-Tribromophenol	68%		39-149%
4165-60-0	Nitrobenzene-d5	73%		32-128%
321-60-8	2-Fluorobiphenyl	68%		35-119%
1718-51-0	Terphenyl-d14	39%		10-126%





MDL = Method Detection Limit

RL = Reporting Limit

E = Indicates value exceeds calibration range

J = Indicates an estimated value

B = Indicates analyte found in associated method blank



SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID: Lab Sample ID:

RA7-GWS JC21103-1

Analyzed

06/01/16

Matrix:

AQ - Ground Water

Date Sampled: Date Received: 05/25/16

05/27/16

Method: Project:

SW846 8270D BY SIM SW846 3510C BMSMC, Building 5 Area, PR

Percent Solids:

n/a

Run #1

File ID DF

By LK Prep Date 05/27/16

Prep Batch OP94327A

Analytical Batch E4P880

Run #2

Initial Volume Final Volume

Run #1 Run #2 920 ml

4P16538A.D

1.0 ml

1

CAS No. Compound Result

RL **MDL** Units

Q

91-20-3

4165-60-0

1718-51-0

321-60-8

Naphthalene

2-Fluorobiphenyl

Terphenyl-d14

ND

0.11

0.032 ug/l

CAS No.

Surrogate Recoveries Nitrobenzene-d5

Run#1

Run# 2 Limits

71% 44% 45% 24-125% 19-127% 10-119%



ND = Not detected

MDL = Method Detection Limit

RL = Reporting Limit

E = Indicates value exceeds calibration range

J = Indicates an estimated value

B = Indicates analyte found in associated method blank



By

XPL

n/a

Page 1 of 1

Client Sample ID: **RA7-GWS** Lab Sample ID: JC21103-1

Matrix: AQ - Ground Water Method: SW846-8015C (DAI)

File ID

GH105253.D

Project: BMSMC, Building 5 Area, PR

DF

1

Date Sampled: 05/25/16

n/a

Q

Date Received: 05/27/16

Percent Solids: n/a

Prep Date **Analytical Batch** Prep Batch

GGH5303

Run #1 Run #2

Low Molecular Alcohol List

CAS No.	Compound	Result	RL	MDL	Units
64-17-5	Ethanol	ND	100	55	ug/l
78-83-1	Isobutyl Alcohol	ND	100	36	ug/l
67-63-0	Isopropyl Alcohol	ND	100	68	ug/l
71-23-8	n-Propyl Alcohol	ND	100	43	ug/l
71-36-3	n-Butyl Alcohol	ND	100	87	ug/l
78-92-2	sec-Butyl Alcohol	ND	100	66	ug/l
67-56-1	Methanol	ND	200	71	ug/l
CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Lim	its
111-27-3	Hexanol	71%		56-1	45%
111-27-3	Hexanol	77%		56-1	45%

Analyzed

05/27/16



ND = Not detected

MDL = Method Detection Limit

RL = Reporting Limit

E = Indicates value exceeds calibration range

J = Indicates an estimated value

B = Indicates analyte found in associated method blank



By

KD

Page 1 of 1

Client Sample ID: RA7-GWS Lab Sample ID:

IC21103-1 AQ - Ground Water Date Sampled: 05/25/16

Matrix: Method:

SW846 8081B SW846 3510C

File ID

4G68800.D

Date Received: 05/27/16

Project:

BMSMC, Building 5 Area, PR

Percent Solids:

Run #1

DF i

Analyzed 06/01/16

Prep Date 05/31/16

Prep Batch OP94371

Q

Analytical Batch G4G1805

Run #2

Initial Volume Final Volume 900 ml

Run #1

10.0 ml

Run #2

Pesticide TCL List

CAS No.	Compound	Result	RL	MDL	Units
309-00-2	Aldrin	ND	0.011	0.0067	ug/l
319-84-6	alpha-BHC	ND	0.011	0.0067	ug/l
319-85-7	beta-BHC	ND	0.011	0.0063	ug/l
319-86-8	delta-BHC	ND	0.011	0.0051	ug/l
58-89-9	gamma-BHC (Lindane)	ND	0.011	0.0031	ug/l
5103-71-9	alpha-Chlordane	ND	0.011	0.0051	ug/l
5103-74-2	gamma-Chlordane	ND	0.011	0.0051	ug/l
60-57-1	Dieldrin	ND	0.011	0.0040	ug/l
72-54-8	4,4'-DDD	ND	0.011	0.0042	ug/l
72-55-9	4,4'-DDE	ND	0.011	0.0068	ug/l
50-29-3	4,4'-DDT	ND	0.011	0.0055	ug/l
72-20-8	Endrin	ND	0.011	0.0056	ug/l
1031-07-8	Endosulfan sulfate	ND	0.011	0.0058	ug/l
7421-93-4	Endrin aldehyde	ND	0.011	0.0057	ug/l
53494-70-5	Endrin ketone	ND	0.011	0.0056	ug/l
959-98-8	Endosulfan-I	ND	0.011	0.0055	ug/l
33213-65-9	Endosulfan-II	ND	0.011	0.0048	ug/l
76-44-8	Heptachlor	ND	0.011	0.0042	ug/l
1024-57-3	Heptachlor epoxide	ND	0.011	0.0073	ug/l
72-43-5	Methoxychlor	ND	0.022	0.0063	ug/l
8001-35-2	Toxaphene	ND	0.28	0.20	ug/l
CAS No.	Surrogate Recoveries	Run#1	Run# 2	Limi	ts
877-09-8	Tetrachloro-m-xylene	76%		26-13	2%
877-09-8	Tetrachloro-m-xylene	84%		26-13	2%
2051-24-3	Decachlorobiphenyl	40%		10-11	8%
2051-24-3	Decachlorobiphenyl	37%		10-11	8%



ND = Not detected

MDL = Method Detection Limit

RL = Reporting Limit

E = Indicates value exceeds calibration range

J = Indicates an estimated value

B = Indicates analyte found in associated method blank



Page 1 of 3

Client Sample ID:	EB05252016
Lab Sample ID:	TC21103-2

Matrix: Method: AQ - Equipment Blank

SW846 8270D SW846 3510C

Date Sampled: 05/25/16 Date Received: 05/27/16

Q

Percent Solids:

Project: BMSMC, Building 5 Area, PR

	File ID	DF	Analyzed	Ву	Prep Date	Prep Batch	Analytical Batch
Run #1	Z111123.D	1	05/28/16	IJ	05/27/16	OP94327	EZ5550
ln 42							

Kun #2

Initial Volume Final Volume 1000 ml $1.0 \, ml$

Run #1 Run #2

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units	
95-57-8	2-Chlorophenol	ND	5.0	0.82	ug/l	
59-50-7	4-Chloro-3-methyl phenol	ND	5.0	0.89	ug/l	
120-83-2	2,4-Dichlorophenol	ND	2.0	1.3	ug/l	
105-67-9	2,4-Dimethylphenol	ND	5.0	2.4	ug/l	
51-28-5	2,4-Dinitrophenol	ND	10	1.6	ug/l	
534-52-1	4,6-Dinitro-o-cresol	ND	5.0	1.3	ug/l	
95-48-7	2-Methylphenol	ND	2.0	0.89	ug/l	
	3&4-Methylphenol	ND	2.0	0.88	ug/l	
88-75-5	2-Nitrophenol	ND	5.0	0.96	ug/l	
100-02-7	4-Nitrophenol	ND	10	1.2	ug/l	
87-86-5	Pentachlorophenol	ND	5.0	1.4	ug/l	
108-95-2	Phenol	ND	2.0	0.39	ug/l	
58-90-2	2,3,4,6-Tetrachlorophenol	ND	5.0	1.5	ug/l	
95-95-4	2,4,5-Trichlorophenol	ND	5.0	1.3	ug/l	
88-06-2	2,4,6-Trichlorophenol	ND	5.0	0.92	ug/l	
83-32-9	Acenaphthene	ND	1:0	0.19	ug/l	
208-96-8	Acenaphthylene	ND	1.0	0.14	ug/l	
98-86-2	Acetophenone	ND	2.0	0.21	ug/l	
120-12-7	Anthracene	ND	1.0	0.21	ug/l	
1912-24-9	Atrazine	ND	2.0	0.45	ug/l	
100-52-7	Benzaldehyde	ND	5.0	0.29	ug/l	
56-55-3	Benzo(a)anthracene	ND	1.0	0.20	ug/l	
50-32-8	Benzo(a)pyrene	ND	1.0	0.21	ug/l	
205-99-2	Benzo(b)fluoranthene	ND	1.0	0.21	ug/I	
191-24-2	Benzo(g,h,i)perylene	ND	1.0	0.34	ug/l	
207-08-9	Benzo(k)fluoranthene	ND	1.0	0.21	ug/l	
101-55-3	4-Bromophenyl phenyl ether	ND	2.0	0.40	ug/l	
85-68-7	Butyl benzyl phthalate	ND	2.0	0.46	ug/l	
92-52-4	1,1'-Biphenyl	ND	1.0	0.21	ug/l	
91-58-7	2-Chloronaphthalene	ND	2.0	0.24	ug/l	
106-47-8	4-Chloroaniline	ND	5.0	0.34	ug/l	
86-74-8	Carbazole	ND	1.0	0.23	ug/l	
					_	



ND = Not detected

MDL = Method Detection Limit

RL = Reporting Limit

E = Indicates value exceeds calibration range

J = Indicates an estimated value

B = Indicates analyte found in associated method blank

Client Sample ID:	EB05252016
Lab Sample ID:	JC21103-2

Matrix: AQ - Equipment Blank

Method: SW846 8270D SW846 3510C Project: BMSMC, Building 5 Area, PR

05/25/16 Date Sampled: Date Received: 05/27/16

Percent Solids: n/a

Q

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units
105-60-2	Caprolactam	ND	2.0	0.65	ug/l
218-01-9	Chrysene	ND	1.0	0.18	ug/l
111-91-1	bis(2-Chloroethoxy)methane	ND	2.0	0.28	ug/l
111-44-4	bis(2-Chloroethyl)ether	ND	2.0	0.25	ug/l
108-60-1	bis(2-Chloroisopropyl)ether	ND	2.0	0.40	ug/l
7005-72-3	4-Chlorophenyl phenyl ether	ND	2.0	0.37	ug/l
121-14-2	2,4-Dinitrotoluene	ND	1.0	0.55	ug/I
606-20-2	2,6-Dinitrotoluene	ND	1.0	0.48	ug/l
91-94-1	3,3'-Dichlorobenzidine	ND	2.0	0.51	ug/l
53-70-3	Dibenzo(a,h)anthracene	ND	1.0	0.33	ug/l
132-64-9	Dibenzofuran	ND	5.0	0.22	ug/l
84-74-2	Di-n-butyl pltthalate	ND	2.0	0.50	ug/l
117-84-0	Di-n-octyl phthalate	ND	2.0	0.23	ug/l
84-66-2	Diethyl phthalate	ND	2.0	0.26	ug/l
131-11-3	Dimethyl phthalate	ND	2.0	0.22	ug/l
117-81-7	bis(2-Ethylhexyl)phthalate	ND	2.0	1.7	ug/l
206-44-0	Fluoranthene	ND	1.0	0.17	ug/l
86-73-7	Fluorene	ND	1.0	0.17	ug/l
118-74-1	Hexachlorobenzene	ND	1.0	0.33	ug/l
87-68-3	Hexachlorobutadiene	ND	1.0	0.49	ug/l
77-47-4	Hexachlorocyclopentadiene	ND	10	2.8	ug/l
67-72-1	Hexachloroethane	ND	2.0	0.39	ug/l
193-39-5	Indeno(1,2,3-cd)pyrene	ND	1.0	0.33	ug/l
78-59-1	Isophorone	ND	2.0	0.28	ug/l
90-12-0	1-Methylnaphthalene	ND	1.0	0.26	ug/l
91-57-6	2-Methylnaphthalene	ND	1.0	0.21	ug/l
88-74-4	2-Nitroaniline	ND	5.0	0.28	ug/l
99-09-2	3-Nitroaniline	ND	5.0	0.39	ug/l
100-01-6	4-Nitroaniline	ND	5.0	0.44	ug/l
98-95-3	Nitrobenzene	ND	2.0	0.64	ug/l
621-64-7	N-Nitroso-di-n-propylamine	ND	2.0	0.48	ug/l
86-30-6	N-Nitrosodiphenylamine	ND	5.0	0.22	ug/l
85-01-8	Phenanthrene	ND	1.0	0.18	ug/l
129-00-0	Pyrene	ND	1.0	0.22	ug/l
95-94-3	1,2,4,5-Tetrachlorobenzene	ND	2.0	0.37	ug/l
					_
CAS No.	Surrogate Recoveries	Run#1	Run# 2	Limi	ts
367-12-4	2-Fluorophenol	37%		14-88	3%
4165-62-2	Phenol-d5	27%		10-11	10%



ND = Not detected MDL = Method Detection Limit

RL = Reporting Limit

E = Indicates value exceeds calibration range

J = Indicates an estimated value

 $B = Indicates \ analyte \ found \ in \ associated \ method \ blank$

N = Indicates presumptive evidence of a compound

Client Sample ID: Lab Sample ID: Matrix:

EB05252016 JC21103-2

AQ - Equipment Blank

Date Sampled: Date Received: 05/25/16 05/27/16

Method: Project:

SW846 8270D SW846 3510C BMSMC, Building 5 Area, PR

Percent Solids: n/a

ABN TCL Special List

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
118-79-6	2,4,6-Tribromophenol	65%		39-149%
4165-60-0	Nitrobenzene-d5	76%		32-128%
321-60-8	2-Fluorobiphenyl	66%		35-119%
1718-51-0	Terphenyl-d14	68%		10-126%



ND = Not detected

MDL = Method Detection Limit

RL = Reporting Limit

E = Indicates value exceeds calibration range

J = Indicates an estimated value

B = Indicates analyte found in associated method blank

Page 1 of 1

Analytical Batch

E4P880

Client Sample ID:	EB05252016
Lab Sample ID:	JC21103-2

Matrix:

SGS Accutest

AQ - Equipment Blank

Date Sampled: 05/25/16 Date Received: 05/27/16

Method: Project:

SW846 8270D BY SIM SW846 3510C

Percent Solids: n/a

BMSMC, Building 5 Area, PR

File ID DF Analyzed By Prep Date Prep Batch 4P16539.D 1 06/01/16 LK 05/27/16 OP94327A

Run #1 Run #2

l		Initial Volume	Final Volume
F	lun #1	1000 ml	1.0 ml

Run #2

CAS No.	Compound	Result	RL	MDL	Units	Q
91-20-3	Naphthalene	ND	0.10	0.029	ug/l	
123-91-1	1,4-Dioxane	ND	0 ₀ 10	0.049	ug/l	

CAS No.	Surrogate Recoveries	Run#1	Run# 2	Limits
4165-60-0	Nitrobenzene-d5	82%		24-125%
321-60-8	2-Fluorobiphenyl	48%		19-127%
1718-51-0	Terphenyl-d14	82%		10-119%



ND = Not detected

MDL = Method Detection Limit

RL = Reporting Limit

E = Indicates value exceeds calibration range

J = Indicates an estimated value

B = Indicates analyte found in associated method blank

Page 1 of 1

Client Sample ID: EB05252016 Lab Sample ID: JC21103-2

Matrix: AQ - Equipment Blank Method: SW846-8015C (DAI)

Project: BMSMC, Building 5 Area, PR

05/25/16 Date Sampled: Date Received: 05/27/16

Q

Percent Solids: n/a

Ву File ID DF Analyzed Prep Date Prep Batch **Analytical Batch** Run #1 GH105256.D 05/27/16 XPL **GGH5303** 1 n/a n/a Run #2

Low Molecular Alcohol List

CAS No.	Compound	Result	RL	MDL	Units
64-17-5	Ethanol	ND	100	55	ug/l
78-83-1	Isobutyl Alcohol	ND	100	36	ug/l
67-63-0	Isopropyl Alcohol	ND	100	68	ug/l
71-23-8	n-Propyl Alcohol	ND	100	43	ug/l
71-36-3	n-Butyl Alcohol	ND	100	87	ug/l
78-92-2	sec-Butyl Alcohol	ND	100	66	ug/l
67-56-1	Methanol	ND	200	71	ug/l
CAS No.	Surrogate Recoveries	Run#1	Run# 2	Lim	its
111-27-3	Hexanol	77%		56-1	45%
111-27-3	Hexanol	81%		56-1	45%



ND = Not detected

MDL = Method Detection Limit

RL = Reporting Limit

E = Indicates value exceeds calibration range

J = Indicates an estimated value

B = Indicates analyte found in associated method blank



Page 1 of 1

Client Sample ID: EB05252016 Lab Sample ID: JC21103-2

File ID

4G68801.D

Matrix: Method: AQ - Equipment Blank

DF

1

SW846 8081B SW846 3510C BMSMC, Building 5 Area, PR

Analyzed

06/01/16

Date Sampled: 05/25/16

Date Received: 05/27/16 Percent Solids: n/a

By Prep Date Prep Batch Analytical Batch
KD 05/31/16 OP94371 G4G1805

Q

Run #1 Run #2

Project:

Initial Volume
1000 ml
Final Volume
10.0 ml

Run #1 Run #2

Pesticide TCL List

CAS No.	Compound	Result	RL	MDL	Units
309-00-2	Aldrin	ND	0.010	0.0060	ug/l
319-84-6	alpha-BHC	ND	0.010	0.0060	ug/l
319-85-7	beta-BHC	ND	0.010	0.0057	ug/l
319-86-8	delta-BHC *	ND	0.010	0.0046	ug/l
58-89-9	gamma-BHC (Lindane)	ND	0.010	0.0028	ug/l
5103-71-9	alpha-Chlordane	ND	0.010	0.0046	ug/l
5103-74-2	gamma-Chlordane	ND	0.010	0.0046	ug/l
60-57-1	Dieldrin	ND	0.010	0.0036	ug/l
72-54-8	4,4'-DDD	ND	0.010	0.0038	ug/l
72-55-9	4,4'-DDE	ND	0.010	0.0062	ug/l
50-29-3	4,4'-DDT	ND	0.010	0.0050	ug/l
72-20-8	Endrin	ND	0.010	0.0050	ug/l
1031-07-8	Endosulfan sulfate	ND	0.010	0.0053	ug/l
7421-93-4	Endrin aldehyde	ND	0.010	0.0051	ug/l
53494-70-5	Endrin ketone	ND	0.010	0.0051	ug/l
959-98-8	Endosulfan-I	ND	0.010	0.0050	ug/l
33213-65-9	Endosulfan-II	ND	0.010	0.0043	ug/l
76-44-8	Heptachlor	ND	0.010	0.0038	ug/l
1024-57-3	Heptachlor epoxide	ND	0.010	0.0065	ug/l
72-43-5	Methoxychlor	ND	0.020	0.0057	ug/l
8001-35-2	Toxaphene	ND	0.25	0.18	ug/l
CAS No.	Surrogate Recoveries	Run#1	Run# 2	Limi	ts
877-09-8	Tetrachloro-m-xylene	62%		26-13	32%
877-09-8	Tetrachloro-m-xylene	59%		26-13	32%
2051-24-3	Decachlorobiphenyl	29%		10-11	8%
2051-24-3	Decachlorobiphenyl	27%		10-11	8%

tack Infante

Veinde/
10 1888

ND = Not detected

MDL = Method Detection Limit

RL = Reporting Limit

E = Indicates value exceeds calibration range

J = Indicates an estimated value

B = Indicates analyte found in associated method blank



SGS Accutest

Report of Analysis

Page 1 of 3

Client Sample ID: SB103 (2.5-3.5) Lab Sample ID:

JC21103-3 SO - Soil

SW846 8270D SW846 3546

Date Sampled: 05/26/16 Date Received:

05/27/16

Method: Project:

Matrix:

BMSMC, Building 5 Area, PR

Percent Solids: 90.9

Ву File ID DF Analyzed Prep Date Prep Batch **Analytical Batch** Run #1 3E84065.D 1 06/10/16 AN 06/01/16 OP94376 E3E3683

Run #2

Final Volume **Initial Weight** 30.1 g

Run #1

 $1.0 \, ml$

Run #2

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units	Q
95-57-8	2-Chlorophenol	ND	73	18	սց/kg	
59-50-7	4-Chloro-3-methyl phenol	ND	180	22	ug/kg	
120-83-2	2,4-Dichlorophenol	ND	180	31	ug/kg	
105-67-9	2,4-Dimethylphenol	ND	180	65	ug/kg	
51-28-5	2,4-Dinitrophenol	ND	180	140	ug/kg	
534-52-1	4,6-Dinitro-o-cresol	ND	180	39	ug/kg	
95-48-7	2-Methylphenol	ND	73	23	ug/kg	
	3&4-Methylphenol	ND	73	30	ug/kg	
88-75-5	2-Nitrophenol	ND	180	24	ug/kg	
100-02-7	4-Nitrophenol	ND	370	98	ug/kg	
87-86-5	Pentachlorophenol	ND	180	34	ug/kg	
108-95-2	Phenol	ND	73	19	ug/kg	
58-90-2	2,3,4,6-Tetrachlorophenol	ND	180	24	ug/kg	
95-95-4	2,4,5-Trichlorophenol	ND	180	27	ug/kg	
88-06-2	2,4,6-Trichlorophenol	ND	180	22	ug/kg	
83-32-9	Acenaphthene	ND	37	13	ug/kg	
208-96-8	Acenaphthylene	ND	37	19	ug/kg	
98-86-2	Acetophenone	ND	180	7.9	ug/kg	
120-12-7	Anthracene	ND	37	22	ug/kg	
1912-24-9	Atrazine	ND	73	16	ug/kg	
56-55-3	Benzo(a)anthracene	ND	37	10	ug/kg	
50-32-8	Benzo(a) pyrene	ND	37	17	ug/kg	
205-99-2	Benzo(b)fluoranthene	ND	37	16	ug/kg	
191-24-2	Benzo(g,h,i)perylene	ND	37	18	ug/kg	
207-08-9	Benzo(k)fluoranthene	ND	37	17	ug/kg	
101-55-3	4-Bromophenyl phenyl ether	ND	73	14	ug/kg	
85-68-7	Butyl benzyl phthalate	ND	73	8.9	ug/kg	
92-52-4	1,1'-Biphenyl	ND	73	5.0	ug/kg	
100-52-7	Benzaldehyde	ND	180	9.1	ug/kg	
91-58-7	2-Chloronaphthalene	ND	73	8.7	ug/kg	
106-47-8	4-Chloroaniline	ND	180	13	ug/kg	
86-74-8	Carbazole	ND	73	5.3	ug/kg	



ND = Not detected

MDL = Method Detection Limit

RL = Reporting Limit

E = Indicates value exceeds calibration range

J = Indicates an estimated value

B = Indicates analyte found in associated method blank



	Sample ID:	SI
Tab C	male ID:	17

D: SB103 (2.5-3.5)

Lab Sample

JC21103-3

SO - Soil

Date Sampled: 05/26/16 Date Received: 05/27/16 Percent Solids: 90.9

Matrix: Method:

SW846 8270D SW846 3546

Project: BMSMC, Building 5 Area, PR

ABN TCL Special List

	-					
CAS No.	Compound	Result	RL	MDL	Units	Q
105-60-2	Caprolactam	ND	73	14	ug/kg	
218-01-9	Chrysene	ND	37	12	ug/kg	
111-91-1	bis(2-Chloroethoxy)methane	ND	73	7.8	ug/kg	
111-44-4	bis(2-Chloroethyl)ether	ND	73	16	ug/kg	
108-60-1	bis(2-Chloroisopropyl)ether	ND	73	13	ug/kg	
7005-72-3	4-Chlorophenyl phenyl ether	ND	73	12	ug/kg	
121-14-2	2,4-Dinitrotoluene	ND	37	11	ug/kg	
606-20-2	2,6-Dinitrotoluene	ND	37	18	ug/kg	
91-94-1	3,3'-Dichlorobenzidine	ND	73	30	ug/kg	
53-70-3	Dibenzo(a,h)anthracene	ND	37	16	ug/kg	
132-64-9	Dibenzofuran	ND	73	15	ug/kg	
84-74-2	Di-n-butyl phthalate	ND	73	6.0	ug/kg	
117-84-0	Di-n-octyl phthalate	ND	73	9.1	ug/kg	
84-66-2	Diethyl phthalate	ND	73	7.8	ug/kg	
131-11-3	Dimethyl phthalate	ND	73	6.5	ug/kg	
117-81-7	bis(2-Ethylhexyl)phthalate	ND	73	8.6	ug/kg	
206-44-0	Fluoranthene	ND	37	16	ug/kg	
86-73-7	Fluorene	ND	37	17	ug/kg	
118-74-1	Hexachlorobenzene	ND	73	9.2	ug/kg	
87-68-3	Hexachlorobutadiene	ND	37	15	ug/kg	
77-47-4	Hexachlorocyclopentadiene	ND	370	15	ug/kg	
67-72-1	Hexachloroethane	ND	180	18	ug/kg	
193-39-5	Indeno(1,2,3-cd)pyrene	ND	37	17	ug/kg	
78-59-1	Isophorone	ND	73	7.8	ug/kg	
90-12-0	1-Methylnaphthalene	ND	73	7.2	ug/kg	
91-57-6	2-Methylnaphthalene	ND	73	8.3	ug/kg	
88-74-4	2-Nitroaniline	ND	180	8.6	ug/kg	
99-09-2	3-Nitroaniline	ND	180	9.1	ug/kg	
100-01-6	4-Nitroaniline	ND	180	9.5	ug/kg	
98-95-3	Nitrobenzene	ND	73	14	ug/kg	HAL ASOCIADO OF PARTIES AND AUGUST AND AUGUS
621-64-7	N-Nitroso-di-n-propylamine	ND	73	11	ug/kg	SOCIADO DE
86-30-6	N-Nitrosodiphenylamine	ND	180	13	ug/kg	Age: Call
85-01-8	Phenanthrene	ND	37	12	ug/kg	(3)
129-00-0	Pyrene	ND	37	12	ug/kg	lact Infante
95-94-3	1,2,4,5-Tetrachlorobenzene	ND	180	9.3	ug/kg	\text{Viendez} \text{Sectors}
CAS No.	Surrogate Recoveries	Run#1	Run# 2	Lim	its	CO LICENCINO
367-12-4	2-Fluorophenol	56%		30-1	106%	COFICENCY
4165-62-2	Phenol-d5	61%			106%	_



MDL = Method Detection Limit

RL = Reporting Limit

E = Indicates value exceeds calibration range

J = Indicates an estimated value

B = Indicates analyte found in associated method blank

N = Indicates presumptive evidence of a compound

Client Sample ID: Lab Sample ID:

SB103 (2.5-3.5) JC21103-3

Matrix:

SO - Soil

Method: Project:

SW846 8270D SW846 3546 BMSMC, Building 5 Area, PR Date Sampled: Date Received;

05/26/16 05/27/16

Percent Solids: 90.9

ABN TCL Special List

CAS No.	Surrogate Recoveries	Run#1	Run# 2	Limits
118-79-6	2,4,6-Tribromophenol	58%		24-140%
4165-60-0	Nitrobenzene-d5	66%		26-122%
321-60-8	2-Fluorobiphenyl	70%		36-112%
1718-51-0	Terphenyl-d14	76%		36-132%



ND = Not detected

MDL = Method Detection Limit

RL = Reporting Limit

E = Indicates value exceeds calibration range

J = Indicates an estimated value

B = Indicates analyte found in associated method blank

Page 1 of 1

Client Sample ID: SB103 (2.5-3.5) Lab Sample ID:

JC21103-3

SO - Soil

Date Sampled: Date Received:

05/26/16 05/27/16

Method: Project:

Matrix:

SW846 8270D BY SIM SW846 3546 BMSMC, Building 5 Area, PR

Percent Solids:

O

90.9

File ID DF Analyzed By Prep Date Prep Batch **Analytical Batch** Run #1 4P16659.D 1 06/04/16]] 06/01/16 OP94376A E4P886

RL

Run #2

Initial Weight Final Volume Run #1 30.1 g

Run #2

123-91-1

91-20-3

1.0 ml

CAS No.

Compound Result 1,4-Dioxane ND Naphthalene ND

3.7 3.7 0.45

MDL

0.73 ug/kg ug/kg

Units

CAS No. Surrogate Recoveries Run#1 Run#2 Limits

4165-60-0 Nitrobenzene-d5 76% 15-138% 12-148% 321-60-8 2-Fluorobiphenyl 64% 1718-51-0 Terphenyl-d14 75% 10-157%



ND = Not detected

MDL = Method Detection Limit

RL = Reporting Limit

E = Indicates value exceeds calibration range

J = Indicates an estimated value

B = Indicates analyte found in associated method blank

N = Indicates presumptive evidence of a compound



JC21103

By

XPL

Prep Date

52-141%

52-141%

n/a

Page 1 of 1

Client Sample ID: Lab Sample ID:

SB103 (2.5-3.5)

Matrix:

JC21103-3 SO - Soil

SW846-8015C (DAI)

Method: Project:

BMSMC, Building 5 Area, PR

DF

Date Sampled: Date Received:

n/a

Q

05/26/16 05/27/16

Percent Solids:

Prep Batch

90.9

Analytical Batch

GGH5305

Run #1 Run #2

Initial Weight

GH105277.D

File ID

Run #1 Run #2

111-27-3

111-27-3

Low Molecular Alcohol List

Hexanol

Hexanol

5.0 g

CAS No.	Compound	Result	RL	MDL	Units	(
64-17-5 78-83-1 67-63-0 71-23-8 71-36-3 78-92-2 67-56-1	Ethanol Isobutyl Alcohol Isopropyl Alcohol n-Propyl Alcohol n-Butyl Alcohol sec-Butyl Alcohol Methanol	ND ND ND ND ND ND	110 110 110 110 110 110 220	76 65 63 44 60 59	ug/kg ug/kg ug/kg ug/kg ug/kg ug/kg ug/kg	
CAS No.	Surrogate Recoveries	Run#1	Run# 2	Lim		

91%

107%

Analyzed

05/31/16



ND = Not detected

MDL = Method Detection Limit

RL = Reporting Limit

E = Indicates value exceeds calibration range

J = Indicates an estimated value

B = Indicates analyte found in associated method blank



Ву

AC

06/01/16

Prep Date

05/31/16

Page 1 of 3

Client Sample ID: SB103D (2.5-3.5)

Lab Sample ID: JC21103-4

Matrix: Method: SO - Soil

6P26567.D

SW846 8270D SW846 3546

Date Sampled: Date Received:

05/26/16 05/27/16

Percent Solids: 91.2

Project: BMSMC, Building 5 Area, PR

File ID DF Analyzed

Prep Batch **Analytical Batch** OP94348 E6P1234

Run #1 Run #2

> **Initial Weight Final Volume**

Run #1 31.4 g 1.0 ml

1

Run #2

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units	Q
95-57-8	2-Chlorophenol	ND	70	17	ug/kg	
59-50-7	4-Chloro-3-methyl phenol	ND	170	21	ug/kg	
120-83-2	2,4-Dichlorophenol	ND	170	30	ug/kg	
105-67-9	2,4-Dimethylphenol	ND	170	62	ug/kg	
51-28-5	2,4-Dinitrophenol	ND	170	130	ug/kg	
534-52-1	4,6-Dinitro-o-cresol	ND	170	37	ug/kg	
95-48-7	2-Methylphenol	ND	70	22	ug/kg	
	3&4-Methylphenol	ND	70	29	ug/kg	
88-75-5	2-Nitrophenol	ND	170	23	ug/kg	
100-02-7	4-Nitrophenol	ND	350	93	ug/kg	
87-86-5	Pentachlorophenol	ND	170	33	ug/kg	
108-95-2	Phenol	ND	70	18	ug/kg	
58-90-2	2,3,4,6-Tetrachlorophenol	ND	170	23	ug/kg	
95-95-4	2,4,5-Trichlorophenol	ND	170	26	ug/kg	
88-06-2	2,4,6-Trichlorophenol	ND	170	21	ug/kg	
83-32-9	Acenaphthene	ND	35	12	ug/kg	
208-96-8	Acenaphthylene	ND	35	18	ug/kg	
98-86-2	Acetophenone	ND	170	7.5	ug/kg	
120-12-7	Anthracene	ND	3 5	21	ug/kg	
1912-24-9	Atrazine	ND	70	15	ug/kg	
56-55-3	Benzo(a)anthracene	ND	35	9.9	ug/kg	
50-32-8	Benzo(a)pyrene	ND	35	16	ug/kg	
205-99-2	Benzo(b)fluoranthene	ND	35	15	ug/kg	
191-24-2	Benzo(g,h,i)perylene	ND	35	17	ug/kg	
207-08-9	Benzo(k)fluoranthene	ND	35	16	ug/kg	
101-55-3	4-Bromophenyl phenyl ether	ND	70	13	ug/kg	
85-68-7	Butyl benzyl phthalate	ND	70	8.5	ug/kg	
92-52-4	1,1'-Biphenyl	ND	70	4.8	ug/kg	
100-52-7	Benzaldehyde	ND	170	8.7	ug/kg	
91-58-7	2-Chloronaphthalene	ND	70	8.3	ug/kg	
106-47-8	4-Chloroaniline	ND	170	13	ug/kg	
86-74-8	Carbazole	ND	70	5.1	ug/kg	



ND = Not detected

MDL = Method Detection Limit

RL = Reporting Limit

E = Indicates value exceeds calibration range

J = Indicates an estimated value

B = Indicates analyte found in associated method blank



Client Sample ID:

SB103D (2.5-3.5)

Lab Sample ID:

JC21103-4

Matrix: Method:

Project:

SO - Soil

SW846 8270D SW846 3546 BMSMC, Building 5 Area, PR Date Sampled: Date Received:

Q

05/26/16 05/27/16

Percent Solids: 91.2

ABN TCL Special List

711A(10D)	openial Disc					
CAS No.	Compound	Result	RL	MDL	Units	
105-60-2	Caprolactam	ND	70	14	ug/kg	
218-01-9	Chrysene	ND	35	11	ug/kg	
111-91-1	bis(2-Chloroethoxy)methane	ND	70	7.5	ug/kg	
111-44-4	bis(2-Chloroethyl)ether	ND	70	15	ug/kg	
108-60-1	bis(2-Chloroisopropyl)ether	ND	70	13	ug/kg	
7005-72-3	4-Chlorophenyl phenyl ether	ND	70	11	ug/kg	
121-14-2	2,4-Dinitrotoluene	ND	35	11	ug/kg	
606-20-2	2,6-Dinitrotoluene	ND	35	18	ug/kg	
91-94-1	3,3'-Dichlorobenzidine	ND	70	29	ug/kg	
53-70-3	Dibenzo(a,h)anthracene	ND	35	15	ug/kg	
132-64-9	Dibenzofuran	ND	70	14	ug/kg	
84-74-2	Di-n-butyl phthalate	ND	70	5.7	ug/kg	
117-84-0	Di-n-octyl phthalate	ND	70	8.7	ug/kg	
84-66-2	Diethyl phthalate	ND	70	7.4	ug/kg	
131-11-3	Dimethyl phthalate	ND	70	6.2	ug/kg	
117-81-7	bis(2-Ethylhexyl)phthalate	ND	70	8.2	ug/kg	
206-44-0	Fluoranthene	ND	35	16	ug/kg	
86-73-7	Fluorene	ND	35	16	ug/kg	
118-74-1	Hexachlorobenzene	ND	70	8.8	ug/kg	
87-68-3	Hexachlorobutadiene	ND	35	14	ug/kg	
77-47-4	Hexachlorocyclopentadiene	ND	350	14	ug/kg	
67-72-1	Hexachloroethane	ND	170	17	ug/kg	
193-39-5	Indeno(1,2,3-cd)pyrene	ND	35	16	ug/kg	
78-59-1	Isophorone	ND	70	7.5	ug/kg	
90-12-0	1-Methylnaphthalene	ND	70	6.8	ug/kg	
91-57-6	2-Methylnaphthalene	ND	70	7.9	ug/kg	
88-74-4	2-Nitroaniline	ND	170	8.2	ug/kg	
99-09-2	3-Nitroaniline	ND	170	8.7	ug/kg	
100-01-6	4-Nitroaniline	ND	170	9.0	ug/kg	
98-95-3	Nitrobenzene	ND	70	13	ug/kg	
621-64-7	N-Nitroso-di-n-propylamine	ND	70	10	ug/kg	
86-30-6	N-Nitrosodiphenylamine	ND	170	13	ug/kg	
85-01-8	Phenanthrene	ND	35	12	ug/kg	
129-00-0	Pyrene	ND	35	11	ug/kg	
95-94-3	1,2,4,5-Tetrachlorobenzene	ND	170	8.9	ug/kg	
CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Lim	its	
367-12-4	2-Fluorophenol	54%		30-1	06%	
4165-62-2	Phenol-d5	61%	30-106%			
NID - Not	ND - Not detected MDI - Method Detection Limit I - Indicates an actim					



ND = Not detected

MDL = Method Detection Limit

RL = Reporting Limit

E = Indicates value exceeds calibration range

J = Indicates an estimated value

B = Indicates analyte found in associated method blank

Client Sample ID: Lab Sample ID:

SB103D (2.5-3.5)

JC21103-4

Matrix: Method:

Project:

SO - Soil

SW846 8270D SW846 3546 BMSMC, Building 5 Area, PR Date Sampled:

05/26/16

Date Received: 05/27/16 Percent Solids: 91.2

ABN TCL Special List

CAS No.	Surrogate Recoveries	Run#1	Run# 2	Limits
118-79-6	2,4,6-Tribromophenol	58%		24-140%
4165-60-0	Nitrobenzene-d5	79%		26-122%
321-60-8	2-Fluorobiphenyl	82%		36-112%
1718-51-0	Terphenyl-d14	88%		36-132%





MDL = Method Detection Limit

RL = Reporting Limit

E = Indicates value exceeds calibration range

J = Indicates an estimated value

B = Indicates analyte found in associated method blank



By

IJ

RL

Run# 2

Page 1 of 1

Client Sample ID: Lab Sample ID:

SB103D (2.5-3.5)

Matrix:

JC21103-4

SO - Soil

06/01/16

05/26/16 Date Sampled:

Date Received:

05/27/16

Method:

SW846 8270D BY SIM SW846 3546

Percent Solids:

Q

91.2

Project:

BMSMC, Building 5 Area, PR

Prep Date

Units

Prep Batch OP94376A

Analytical Batch E4P886

Run #1 Run #2

> Initial Weight **Final Volume** 30.6 g

File ID

4P16660.D

Run #1 Run #2

91-20-3

1.0 ml

DF

1

CAS No.	Compound
123-91-1	1,4-Dioxane a

Naphthalene

ND ND

Result

Analyzed

06/04/16

3.6 0.72 3.6 0.44

MDL

ug/kg ug/kg

CAS No. Surrogate Recoveries

4165-60-0 Nitrobenzene-d5 321-60-8 2-Fluorobiphenyl 1718-51-0 Terphenyl-d14

Run#1 84%

58%

87%

15-138% 12-148% 10-157%

Limits

(a) Not accredited for this compound at the time of analysis, but all method requirements were followed.



ND = Not detected

E = Indicates value exceeds calibration range

J = Indicates an estimated value

B = Indicates analyte found in associated method blank

N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID: SB103D (2.5-3.5) Lab Sample ID: JC21103-4

File ID

Matrix: Method:

Project:

SO - Soil

SW846-8015C (DAI)

DF

1

Date Sampled: 05/26/16 Date Received: 05/27/16

BMSMC, Building 5 Area, PR

Percent Solids: 91.2

Ву Analyzed **Analytical Batch** Prep Date Prep Batch 05/31/16 XPL GGH5305 n/a n/a

Run #1 Run #2

Initial Weight

GH105280.D

Run #1

Run #2

111-27-3

Low Molecular Alcohol List

Hexanol

5.1 g

CAS No.	Compound	Result	RL	MDL	Units	Q
64-17-5	Ethanol	ND	110	74	ug/kg	
78-83-1	Isobutyl Alcohol	ND	110	63	ug/kg	
67-63-0	Isopropyl Alcohol	ND	110	61	ug/kg	
71-23-8	n-Propyl Alcohol	ND	110	43	ug/kg	
71-36-3	n-Butyl Alcohol	ND	110	58	ug/kg	
78-92-2	sec-Butyl Alcohol	ND	110	57	ug/kg	
67-56-1	Methanol	ND	210	51	ug/kg	
CAS No.	Surrogate Recoveries	Run#1	Run# 2	Lim	its	
111-27-3	Hexanol	94%		52-1	41%	

99%



ND = Not detected

MDL = Method Detection Limit

RL = Reporting Limit

E = Indicates value exceeds calibration range

J = Indicates an estimated value

52-141%

B = Indicates analyte found in associated method blank

SB103 (6-7) JC21103-5

Matrix:

SO - Soil

SW846 8270D SW846 3546

DF

1

Date Sampled: 05/26/16 Date Received: 05/27/16

Method: Project:

BMSMC, Building 5 Area, PR

Percent Solids: 85.6

Run #1

File ID 3P54297.D Analyzed 06/13/16

Ву RL Prep Date 06/01/16

Prep Batch OP94376

Analytical Batch E3P2473

Run #2

Initial Weight

Final Volume

Run #1

30.4 g

1.0 ml

Run #2

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units	Q
95-57-8	2-Chlorophenol	ND	77	19	ug/kg	
59-50-7	4-Chloro-3-methyl phenol	ND	190	24	ug/kg	
120-83-2	2,4-Dichlorophenol	ND	190	33	ug/kg	
105-67-9	2,4-Dimethylphenol	ND	190	68	ug/kg	
51-28-5	2,4-Dinitrophenol	ND	190	140	ug/kg	
534-52-1	4,6-Dinitro-o-cresol	ND	190	41	ug/kg	
95-48-7	2-Methylphenol	ND	77	25	ug/kg	
	3&4-Methylphenol	ND	77	32	ug/kg	
88-75-5	2-Nitrophenol	ND	190	25	ug/kg	
100-02-7	4-Nitrophenol	ND	380	100	ug/kg	
87-86-5	Pentachlorophenol	ND	190	36	ug/kg	
108-95-2	Phenoi	ND	77	20	ug/kg	
58-90-2	2,3,4,6-Tetrachlorophenol	ND	190	25	ug/kg	
95-95-4	2,4,5-Trichlorophenol	ND	190	29	ug/kg	
88-06-2	2,4,6-Trichlorophenol	ND	190	23	ug/kg	
83-32-9	Acenaphthene	ND	38	13	ug/kg	
208-96-8	Acenaphthylene	ND	38	20	ug/kg	
98-86-2	Acetophenone	ND	190	8.3	ug/kg	
120-12-7	Anthracene	56.8	38	24	ug/kg	
1912-24-9	Atrazine	ND	77	16	ug/kg	
56-55-3	Benzo(a) anthracene	76.3	38	11	ug/kg	
50-32-8	Benzo(a)pyrene	48.7	38	17	ug/kg	
205-99-2	Benzo(b)fluoranthene	66.4	38	17	ug/kg	
191-24-2	Benzo(g,h,i)perylene	19.8	38	19	ug/kg	J
207-08-9	Benzo(k)fluoranthene	24.9	38	18	ug/kg	j
101-55-3	4-Bromophenyl phenyl ether	ND	77	15	ug/kg	-
85-68-7	Butyl benzyl phthalate	ND	77	9.4	ug/kg	
92-52-4	1,1'-Biphenyl	ND	77	5.3	ug/kg	
100-52-7	Benzaldehyde	ND	190	9.5	ug/kg	
91-58-7	2-Chloronaphthalene	ND	77	9.1	ug/kg	
106-47-8	4-Chloroaniline	ND	190	14	ug/kg	
					0 0	

taci infante \téndez

ND = Not detected

86-74-8

MDL = Method Detection Limit

ND

77

5.6

RL = Reporting Limit

E = Indicates value exceeds calibration range

Carbazole

J = Indicates an estimated value

ug/kg

B = Indicates analyte found in associated method blank



Client Sample ID: Lab Sample ID:

SB103 (6-7) JC21103-5

Matrix:

SO - Soil

Method:

SW846 8270D SW846 3546

Project:

BMSMC, Building 5 Area, PR

Date Sampled: 05/26/16

Date Received: 05/27/16

Q

J

j

Percent Solids: 85.6

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units
105-60-2	Caprolactam	ND	77	15	ug/kg
218-01-9	Chrysene	25.9	38	12	ug/kg
111-91-1	bis(2-Chloroethoxy)methane	ND	77	8.2	ug/kg
111-44-4	bis(2-Chloroethyl)ether	ND	77	17	ug/kg
108-60-1	bis(2-Chloroisopropyl)ether	ND	77	14	ug/kg
7005-72-3	4-Chlorophenyl phenyl ether	ND	77	12	ug/kg
121-14-2	2,4-Dinitrotoluene	ND	38	12	ug/kg
606-20-2	2,6-Dinitrotoluene	ND	38	19	ug/kg
91-94-1	3,3'-Dichlorobenzidine	ND	77	32	ug/kg
53-70-3	Dibenzo(a,h)anthracene	ND	38	17	ug/kg
132-64-9	Dibenzofuran	ND	77	16	ug/kg
84-74-2	Di-n-butyl phthalate	ND	77	6.3	ug/kg
117-84-0	Di-n-octyl phthalate	NĐ	77	9.6	ug/kg
84-66-2	Diethyl phthalate	ND	77	8.2	ug/kg
131-11-3	Dimethyl phthalate	ND	77	6.8	ug/kg
117-81-7	bis(2-Ethylhexyl)phthalate	ND	77	9.0	ug/kg
206-44-0	Fluoranthene	501	38	17	ug/kg
86-73-7	Fluorene	ND	38	18	ug/kg
118-74-1	Hexachlorobenzene	ND	77	9.7	ug/kg
87-68-3	Hexachlorobutadiene	ND	38	15	ug/kg
77-47-4	Hexachlorocyclopentadiene	ND	380	15	ug/kg
67-72-1	Hexachloroethane	ND	190	19	ug/kg
193-39-5	Indeno(1,2,3-cd)pyrene	21.2	38	18	ug/kg
78-59-1	Isophorone	ND	77	8.2	ug/kg
90-12-0	1-Methylnaphthalene	ND	77	7.5	ug/kg
91-57-6	2-Methylnaphthalene	ND	77	8.7	ug/kg
88-74-4	2-Nitroaniline	ND	190	9.1	ug/kg
99-09-2	3-Nitroaniline	ND	190	9.6	ug/kg
100-01-6	4-Nitroaniline	ND	190	10	ug/kg
98-95-3	Nitrobenzene	ND	77	15	ug/kg
621-64-7	N-Nitroso-di-n-propylamine	ND	77	11	ug/kg
86-30-6	N-Nitrosodiphenylamine	ND	190	14	ug/kg
85-01-8	Phenanthrene	ND	38	13	ug/kg
129-00-0	Pyrene	321	38	12	ug/kg
95-94-3	1,2,4,5-Tetrachlorobenzene	ND	190	9.8	ug/kg
CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limit	s



ND = Not detected

367-12-4

4165-62-2

MDL = Method Detection Limit

86%

81%

RL = Reporting Limit

E = Indicates value exceeds calibration range

2-Fluorophenol

Phenol-d5

J = Indicates an estimated value

30-106%

30-106%

B = Indicates analyte found in associated method blank





Client Sample ID: SB103 (6-7) Lab Sample ID:

JC21103-5

Matrix:

SO - Soil

SW846 8270D SW846 3546

Date Sampled: Date Received:

05/26/16

Percent Solids: 85.6

05/27/16


Method: Project:

BMSMC, Building 5 Area, PR

ABN TCL Special List

CAS No.	Surrogate Recoveries	Run#1	Run# 2	Limits
118-79-6	2,4,6-Tribromophenol	89%		24-140%
4165-60-0	Nitrobenzene-d5	113%		26-122%
321-60-8	2-Fluorobiphenyl	85%		36-112%
1718-51-0	Terphenyl-d14	86%		36-132%



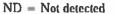

MDL = Method Detection Limit

RL = Reporting Limit

E = Indicates value exceeds calibration range

J = Indicates an estimated value

B = Indicates analyte found in associated method blank





SGS Accutest

Report of Analysis

By

IJ

Page 1 of 1

Client Sample ID: Lab Sample ID:

SB103 (6-7) JC21103-5

Matrix: Method: SO - Soil

SW846 8270D BY SIM SW846 3546

Analyzed

06/04/16

Date Sampled: Date Received:

05/26/16 05/27/16

Percent Solids:

OP94376A

85.6

Project:

BMSMC, Building 5 Area, PR

Units

Prep Date

06/01/16

MDL

0.77

Prep Batch

Q

Analytical Batch E4P886

Run #1 Run #2

91-20-3

Initial Weight Run #1

Final Volume

DF

1

30.4 g

File ID

4P16661.D

 $1.0 \, \mathrm{ml}$

Run #2 CAS No. RLCompound Result 123-91-1 1,4-Dioxane a ND 3.8

CAS No. Surrogate Recoveries

Naphthalene

4165-60-0 Nitrobenzene-d5

321-60-8 2-Fluorobiphenyl 1718-51-0 Terphenyl-d14

ug/kg ND 3.8 0.47 ug/kg

Run#1 Run#2 Limits 70%

15-138% 12-148% 10-157%

(a) Not accredited for this compound at the time of analysis, but all method requirements were followed.

51%

72%



ND = Not detected

MDL = Method Detection Limit

RL = Reporting Limit

E = Indicates value exceeds calibration range

J = Indicates an estimated value

B = Indicates analyte found in associated method blank





SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID: SB103 (6-7) Lab Sample ID:

JC21103-5 SO - Soil

Date Sampled: 05/26/16

Matrix:

SW846-8015C (DAI)

Date Received: 05/27/16

Method:

Percent Solids: 85.6

Project:

BMSMC, Building 5 Area, PR

Run #1

File ID DF GH105281.D 1

Ву Analyzed 05/31/16 XPL

Prep Date n/a

Prep Batch n/a

Analytical Batch GGH5305

Run #2

Initial Weight

Run #1 Run #2

111-27-3

5.1 g

Low Molecular Alcohol List

Hexanol

CAS No.	Compound	Result	RL	MDL	Units	Q
64-17-5	Ethanol	ND	110	79	ug/kg	

64-17-5	Ethanol	ND	110	79	ug/kg
78-83-1	Isobutyl Alcohol	ND	110	67	ug/kg
67-63-0	Isopropyl Alcohol	ND	110	65	ug/kg
71-23-8	n-Propyl Alcohol	ND	110	46	ug/kg
71-36-3	n-Butyl Alcohol	ND	110	62	ug/kg
78-92-2	sec-Butyl Alcohol	ND	110	61	ug/kg
67-56-1	Methanol	ND	230	55	ug/kg

0. 00 1	2720 22322 22	.,,	200	00 56
CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
111-27-3	Hexanol	95%		52-141%

100%



ND = Not detected

MDL = Method Detection Limit

RL = Reporting Limit

E = Indicates value exceeds calibration range

J = Indicates an estimated value

52-141%

B = Indicates analyte found in associated method blank



Page 1 of 3

Client Sample ID: Lab Sample ID:

MW-23S (5-6) JC21103-6

Matrix:

Method:

SO - Soil

SW846 8270D SW846 3546 BMSMC, Building 5 Area, PR Date Sampled:

05/26/16

Date Received: 05/27/16

Percent Solids: 82.8

Run #1 Run #2

Project:

File ID 6P26569.D DF 1

Analyzed Ву 06/01/16 AC Prep Date 05/31/16

Prep Batch OP94348

Q

Analytical Batch E6P1234

Initial Weight Run #1 30.1 g

Final Volume 1.0 ml

Run #2

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units
95-57-8	2-Chlorophenol	ND	80	20	ug/kg
59-50-7	4-Chloro-3-methyl phenol	ND	200	25	ug/kg
120-83-2	2,4-Dichlorophenol	ND	200	34	ug/kg
105-67-9	2,4-Dimethylphenol	ND	200	71	ug/kg
51-28-5	2,4-Dinitrophenol	ND	200	150	ug/kg
534-52-1	4,6-Dinitro-o-cresol	ND	200	43	ug/kg
95-48-7	2-Methylphenol	ND	80	26	ug/kg
	3&4-Methylphenol	ND	80	33	ug/kg
88-75-5	2-Nitrophenol	ND	200	27	ug/kg
100-02-7	4-Nitrophenol	ND	400	110	ug/kg
87-86-5	Pentachlorophenol	ND	200	38	ug/kg
108-95-2	Phenol	ND	80	21	ug/kg
58-90-2	2,3,4,6-Tetrachlorophenol	ND	200	27	ug/kg
95-95-4	2,4,5-Trichlorophenol	ND	200	30	ug/kg
88-06-2	2,4,6-Trichlorophenol	ND	200	24	ug/kg
83-32-9	Acenaphthene	ND	40	14	ug/kg
208-96-8	Acenaphthylene	ND	40	20	ug/kg
98-86-2	Acetophenone	ND	200	8.6	ug/kg
120-12-7	Anthracene	ND	40	25	ug/kg
1912-24-9	Atrazine	ND	80	17	ug/kg
56-55-3	Benzo(a)anthracene	ND	40	11	ug/kg
50-32-8	Вепло(а)ругене	ND	40	18	ug/kg
205-99-2	Benzo(b)fluoranthene	ND	40	18	ug/kg
191-24-2	Benzo(g,h,i)perylene	ND	40	20	ug/kg
207-08-9	Benzo(k)fluoranthene	ND	40	19	ug/kg
101-55-3	4-Bromophenyl phenyl ether	ND	80	15	ug/kg
85-68-7	Butyl benzyl phthalate	ND	80	9.8	ug/kg
92-52-4	1,1'-Biphenyl	ND	80	5.5	ug/kg
100-52-7	Benzaldehyde	ND	200	10	ug/kg

BRE ASOCIADO DE fael Infante Méndez 1888

ND = Not detected

91-58-7

106-47-8

86-74-8

MDL = Method Detection Limit

ND

ND

ND

RL = Reporting Limit

E = Indicates value exceeds calibration range

4-Chloroaniline

Carbazole

2-Chloronaphthalene

J = Indicates an estimated value

ug/kg

ug/kg

ug/kg

9.5

14

5.8

80

200

B = Indicates analyte found in associated method blank



Client Sample ID: MW-23S (5-6) Lab Sample ID:

JC21103-6

Matrix: Method:

Project:

SO - Soil

SW846 8270D SW846 3546

BMSMC, Building 5 Area, PR

Date Sampled:

Q

05/26/16 Date Received: 05/27/16

Percent Solids: 82.8

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units
105-60-2	Caprolactam	ND	80	16	ug/kg
218-01-9	Chrysene	ND	40	13	ug/kg
111-91-1	bis(2-Chloroethoxy)methane	ND	80	8.6	ug/kg
111-44-4	bis(2-Chloroethyl)ether	ND	80	17	ug/kg
108-60-1	bis(2-Chloroisopropyl)ether	ND	80	14	ug/kg
7005-72-3	4-Chlorophenyl phenyl ether	ND	80	13	ug/kg
121-14-2	2,4-Dinitrotoluene	ND	40	12	ug/kg
606-20-2	2,6-Dinitrotoluene	ND	40	20	ug/kg
91-94-1	3,3'-Dichlorobenzidine	ND	80	33	ug/kg
53-70-3	Dibenzo(a,h)anthracene	ND	40	18	ug/kg
132-64-9	Dibenzofuran	ND	80	16	ug/kg
84-74-2	Di-n-butyl phthalate	ND	80	6.5	ug/kg
117-84-0	Di-n-octyl phthalate	ND	80	10	ug/kg
84-66-2	Diethyl phthalate	ND	80	8.5	ug/kg
131-11-3	Dimethyl phthalate	ND	80	7.1	ug/kg
117-81-7	bis(2-Ethylhexyl)phthalate	ND	80	9.4	ug/kg
206-44-0	Fluoranthene	ND	40	18	ug/kg
86-73-7	Fluorene	ND	40	18	ug/kg
118-74-1	Hexachlorobenzene	ND	80	10	ug/kg
87-68-3	Hexachlorobutadiene	ND	40	16	ug/kg
77-47-4	Hexachlorocyclopentadiene	ND	400	16	ug/kg
67-72-1	Hexachloroethane	ND	200	20	ug/kg
193-39-5	Indeno(1,2,3-cd)pyrene	ND	40	19	ug/kg
78-59-1	Isophorone	ND	80	8.6	ug/kg
90-12-0	1-Methylnaphthalene	ND	80	7.9	ug/kg
91-57-6	2-Methylnaphthalene	ND	80	9.1	ug/kg
88-74-4	2-Nitroaniline	ND	200	9.5	ug/kg
99-09-2	3-Nitroaniline	ND	200	10	ug/kg
100-01-6	4-Nitroaniline	ND	200	10	ug/kg
98-95-3	Nitrobenzene	ND	80	15	ug/kg
621-64-7	N-Nitroso-di-n-propylamine	ND	80	12	ug/kg
86-30-6	N-Nitrosodiphenylamine	ND	200	15	ug/kg
85-01-8	Phenanthrene	ND	40	13	ug/kg
129-00-0	Pyrene	ND	40	13	ug/kg
95-94-3	1,2,4,5-Tetrachlorobenzene	ND	200	10	ug/kg
CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limi	its



ND = Not detected

367-12-4

4165-62-2

MDL = Method Detection Limit

72%

71%

RL = Reporting Limit

E = Indicates value exceeds calibration range

2-Fluorophenol

Phenol-d5

J = Indicates an estimated value

30-106%

30-106%

B = Indicates analyte found in associated method blank





Client Sample ID: MW-23S (5-6) Lab Sample ID:

JC21103-6 SO - Soil

SW846 8270D SW846 3546

Date Sampled: Date Received: 05/27/16 Percent Solids: 82.8

05/26/16

Method: Project:


Matrix:

BMSMC, Building 5 Area, PR

ABN TCL Special List

CAS No.	Surrogate Recoveries	Run#1	Run# 2	Limits
118-79-6	2,4,6-Tribromophenol	71%		24-140%
4165-60-0	Nitrobenzene-d5	80%		26-122%
321-60-8	2-Fluorobiphenyl	83%		36-112%
1718-51-0	Terphenyl-d14	92%		36-132%



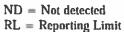

MDL = Method Detection Limit

E = Indicates value exceeds calibration range

J = Indicates an estimated value

B = Indicates analyte found in associated method blank





By

Ш

Prep Date

10-157%

06/01/16

Page 1 of 1

Client Sample ID: Lab Sample ID:

MW-23S (5-6) JC21103-6

SO - Soil

Date Sampled: Date Received:

05/26/16 05/27/16

Matrix: Method:

SW846 8270D BY SIM SW846 3546

Analyzed

06/04/16

Percent Solids: 82.8

Project:

BMSMC, Building 5 Area, PR

DF

1

Prep Batch

OP94376A

Analytical Batch E4P886

Run #1 Run #2

Initial Weight Final Volume

Terphenyl-d14

File ID

4P16662.D

Run #1

30.7 g1.0 ml

Run #2

1718-51-0

CAS No. Compound Result RL MDL Units Q 123-91-1 1,4-Dioxane a ND 3.9 0.79 ug/kg 91-20-3 Naphthalene ND 3.9 0.48 ug/kg CAS No. Surrogate Recoveries Run#1 Run# 2 Limits 4165-60-0 Nitrobenzene-d5 84% 15-138% 321-60-8 2-Fluorobiphenyl 58% 12-148%

(a) Not accredited for this compound at the time of analysis, but all method requirements were followed.

73%



ND = Not detected

MDL = Method Detection Limit

RL = Reporting Limit

E = Indicates value exceeds calibration range

J = Indicates an estimated value

B = Indicates analyte found in associated method blank



SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID: Lab Sample ID:

MW-23S (5-6) JC21103-6

Matrix:

SO - Soil

Method:

SW846-8015C (DAI)

Project:

BMSMC, Building 5 Area, PR

Date Sampled:

05/26/16

Date Received: 05/27/16

Percent Solids: 82.8

File ID DF Analyzed By **Analytical Batch** Prep Date Prep Batch Run #1 GH105282.D 1 05/31/16 XPL GGH5305 n/a n/a

Run #2

Initial Weight

Run #1 5.0 g

Run #2

Low Molecular Alcohol List.

CAS No. Compound Result RL MDL Units Q 64-17-5 Ethanol ND 120 83 ug/kg 78-83-1 Isobutyl Alcohol ND 120 71 ug/kg 67-63-0 Isopropyl Alcohol ug/kg ND 120 69 71-23-8 n-Propyl Alcohol ND 120 49 ug/kg 71-36-3 n-Butyl Alcohol ND 120 66 ug/kg 78-92-2 sec-Butyl Alcohol ND 120 64 ug/kg 67-56-1 Methanol ND 240 58 ug/kg

CAS No. Surrogate Recoveries Run#1 Run#2 Limits 111-27-3 Hexanol 97% 52-141% 111-27-3 Hexanol 103% 52-141%





MDL = Method Detection Limit



RL = Reporting Limit

E = Indicates value exceeds calibration range

^{] =} Indicates an estimated value

B = Indicates analyte found in associated method blank

N = Indicates presumptive evidence of a compound

By

IJ

Page 1 of 3

Client Sample ID: Lab Sample ID:

SB103-GWD JC21103-7

Matrix:

Method:

AQ - Ground Water SW846 8270D SW846 3510C

DF

1

Date Sampled:

05/26/16 Date Received: 05/27/16

Percent Solids: n/a

Prep Date

05/27/16

Project:

BMSMC, Building 5 Area, PR

Analyzed

05/28/16

Prep Batch OP94327

Q

Analytical Batch EZ5550

Run #1 Run #2

Initial Volume Final Volume 900 ml

Run #1 Run #2

File ID

Z111121.D

1.0 ml

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units
95-57-8	2-Chlorophenol	ND	5.6	0.91	ug/l
59-50-7	4-Chloro-3-methyl phenol	ND	5.6	0.99	ug/l
120-83-2	2,4-Dichlorophenol	ND	2.2	1.4	ug/l
105-67-9	2,4-Dimethylphenol	ND	5.6	2.7	ug/l
51-28-5	2,4-Dinitrophenol	ND	11	1.7	ug/l
534-52-1	4,6-Dinitro-o-cresol	ND	5.6	1.4	ug/l
95-48-7	2-Methylphenol	ND	2.2	0.99	ug/l
	3&4-Methylphenol	ND	2.2	0.98	ug/l
88-75-5	2-Nitrophenol	ND	5.6	1.1	ug/l
100-02-7	4-Nitrophenol	ND	11	1.3	ug/l
87-86-5	Pentachlorophenol	ND	5.6	1.5	ug/l
108-95-2	Phenol	ND	2.2	0.44	ug/l
58-90-2	2,3,4,6-Tetrachlorophenol	ND	5.6	1.6	ug/l
95-95-4	2,4,5-Trichlorophenol	ND	5.6	1.5	ug/l
88-06-2	2,4,6-Trichlorophenol	ND	5.6	1.0	ug/l
83-32-9	Acenaphthene	ND	1.1	0.21	ug/l
208-96-8	Acenaphthylene	ND	1.1	0.15	ug/l
98-86-2	Acetophenone	ND	2.2	0.23	ug/l
120-12-7	Anthracene	ND	1.1	0.23	ug/l
1912-24-9	Atrazine	ND	2.2	0.50	ug/l
100-52-7	Benzaldehyde	ND	5.6	0.32	ug/l
56-55-3	Benzo(a)anthracene	ND	1.1	0.23	ug/l
50-32-8	Benzo(a)pyrene	ND	1.1	0.24	ug/l
205-99-2	Benzo(b)fluoranthene	ND	1.1	0.23	ug/l
191-24-2	Benzo(g,h,i)perylene	ND	1.1	0.38	ug/l
207-08-9	Benzo(k)fluoranthene	ND	1.1	0.23	ug/l
101-55-3	4-Bromophenyl phenyl ether	ND	2.2	0.45	ug/l
85-68-7	Butyl benzyl phthalate	ND	2.2	0.51	ug/l
92-52-4	1,1'-Biphenyl	ND	1:1	0.24	ug/l
91-58-7	2-Chloronaphthalene	ND	2.2	0.26	ug/l
106-47-8	4-Chloroaniline	ND	5.6	0.38	ug/l
86-74-8	Carbazole	ND	1.1	0.25	ug/l

taci Infante Méndez

ND = Not detected

MDL = Method Detection Limit

RL = Reporting Limit

E = Indicates value exceeds calibration range

J = Indicates an estimated value

B = Indicates analyte found in associated method blank



Client Sample ID: Lab Sample ID:

SB103-GWD JC21103-7

Matrix:

AQ - Ground Water

Method: Project:

SW846 8270D SW846 3510C BMSMC, Building 5 Area, PR

Date Sampled:

Q

05/26/16 Date Received: 05/27/16

Percent Solids: n/a

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units	
105-60-2	Caprolactam	ND	2.2	0.72	ug/l	
218-01-9	Chrysene	ND	1.1	0.20	ug/l	
111-91-1	bis(2-Chloroethoxy)methane	ND	2.2	0.31	ug/l	
111-44-4	bis(2-Chloroethyl)ether	ND	2.2	0.28	ug/l	
108-60-1	bis(2-Chloroisopropyl)ether	ND	2.2	0.45	ug/l	
7005-72-3	4-Chlorophenyl phenyl ether	ND	2.2	0.41	ug/l	
121-14-2	2,4-Dinitrotoluene	ND	1.1	0.61	ug/l	
606-20-2	2,6-Dinitrotoluene	ND	1.1	0.53	ug/l	
91-94-1	3,3'-Dichlorobenzidine	ND	2.2	0.56	ug/l	
53-70-3	Dibenzo(a,b)anthracene	ND	1.1	0.37	ug/l	
132-64-9	Dibenzofuran	ND	5.6	0.24	ug/l	
84-74-2	Di-n-butyl phthalate	ND	2.2	0.55	ug/l	
117-84-0	Di-n-octyl phthalate	ND	2.2	0.26	ug/l	
84-66-2	Diethyl phthalate	ND	2.2	0.29	ug/l	
131-11-3	Dimethyl phthalate	ND	2.2	0.24	ug/l	
117-81-7	bis(2-Ethylhexyl)phthalate	ND	2.2	1.8	ug/l	
206-44-0	Fluoranthene	ND	1.1	0.19	ug/l	
86-73-7	Fluorene	ND	1.1	0.19	ug/l	
118-74-1	Hexachlorobenzene	ND	1.1	0.36	ug/l	
87-68-3	Hexachlorobutadiene	ND	1.1	0.55	ug/l	
77-47-4	Hexachlorocyclopentadiene	ND	11	3.1	ug/l	
67-72-1	Hexachloroethane	ND	2.2	0.43	ug/l	
193-39-5	Indeno(1,2,3-cd)pyrene	ND	1.1	0.37	ug/l	
78-59-1	Isophorone	ND	2.2	0.31	ug/l	
90-12-0	1-Methylnaphthalene	ND	1.1	0.29	ug/l	
91-57-6	2-Methylnaphthalene	ND	1.1	0.23	ug/l	
88-74-4	2-Nitroaniline	ND	5.6	0.31	ug/l	
99-09-2	3-Nitroaniline	ND	5.6	0.43	ug/l	
100-01-6	4-Nitroaniline	ND	5.6	0.49	ug/l	
98-95-3	Nitrobenzene	ND	2.2	0.71	ug/l	
621-64-7	N-Nitroso-di-n-propylamine	ND	2.2	0.53	ug/l	
86-30-6	N-Nitrosodiphenylamine	ND	5.6	0.25	ug/l	
85-01-8	Phenanthrene	ND	1.1	0.19	ug/l	
129-00-0	Pyrene	ND	1.1	0.24	ug/l	
95-94-3	1,2,4,5-Tetrachlorobenzene	ND	2.2	0.41	ug/l	
CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits		
367-12-4	2-Fluorophenol	36%		14-8	14-88%	
4165-62-2	Phenoi-d5	28%		10-1	10%	



ND = Not detected

MDL = Method Detection Limit

RL = Reporting Limit

J = Indicates an estimated value

B = Indicates analyte found in associated method blank



E = Indicates value exceeds calibration range

Report of Analysis

Page 3 of 3

Client Sample ID: Lab Sample ID:

SB103-GWD JC21103-7

AQ - Ground Water

Matrix: Method: Project:

SW846 8270D SW846 3510C BMSMC, Building 5 Area, PR

Date Sampled: Date Received: 05/27/16

05/26/16

Percent Solids: n/a

ABN TCL Special List

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
118-79-6 4165-60-0	2,4,6-Tribromophenol Nitrobenzene-d5	63% 69%		39-149% 32-128%
321-60-8	2-Fluorobiphenyl	61%		35-119%
1718-51-0	Terphenyl-d14	40%		10-126%





MDL = Method Detection Limit

RL = Reporting Limit

E = Indicates value exceeds calibration range

J = Indicates an estimated value

B = Indicates analyte found in associated method blank





SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID: Lab Sample ID:

SB103-GWD JC21103-7

Matrix:

AQ - Ground Water

DF

SW846 8270D BY SIM SW846 3510C

Date Sampled: Date Received:

05/26/16 05/27/16

Method: Project:

BMSMC, Building 5 Area, PR

Percent Solids: n/a

0

Run #2

Run #1 3M61852.D 1 06/01/16 IJ

By

Analyzed

Prep Date 05/27/16

Prep Batch OP94327A

Analytical Batch E3M2915

Initial Volume **Final Volume** Run #1 900 ml 1.0 ml

Run #2

CAS No. Compound

File ID

91-20-3 Naphthalene 123-91-1 1,4-Dioxane

CAS No. Surrogate Recoveries

4165-60-0 Nitrobenzene-d5 321-60-8 2-Fluorobiphenyl 1718-51-0 Terphenyl-d14

Result RL

Run#1

56%

74%

46%

ND 0.11 0.2980.11

Run# 2

0.033 0.054

MDL

ug/l ug/l

Units

Limits 24-125% 19-127%

10-119%

tacl Infante Ménde/

ND = Not detected

MDL = Method Detection Limit

RL = Reporting Limit

E = Indicates value exceeds calibration range

J = Indicates an estimated value

B = Indicates analyte found in associated method blank



Report of Analysis

Page 1 of 1

Client Sample ID: Lab Sample ID:

SB103-GWD JC21103-7

Matrix: Method:

AQ - Ground Water SW846-8015C (DAI)

Project:

BMSMC, Building 5 Area, PR

Date Sampled: 05/26/16

Date Received: 05/27/16

Percent Solids: n/a

File ID DF Analyzed By Run #1 GH105257.D 1 05/27/16 XPL Run #2		rep Batch Analytical Batch /a GGH5303
--	--	--

Low Molecular Alcohol List

CAS No.	Compound	Result	RL	MDL	Units	Q
64-17-5	Ethanol	ND	100	55	ug/l	
78-83-1	Isobutyl Alcohol	ND	100	36	ug/l	
67-63-0	Isopropyl Alcohol	ND	100	68	ug/l	
71-23-8	n-Propyl Alcohol	ND	100	43	ug/l	
71-36-3	n-Butyl Alcohol	ND	100	87	ug/l	
78-92-2	sec-Butyl Alcohol	ND	100	66	ug/l	
67-56-1	Methanol	ND	200	71	ug/l	
CAS No.	Surrogate Recoveries	Run#1	Run# 2	Lim	its	
111-27-3	Hexanol	88%		56-1	45%	
111-27-3	Hexanol	92%		56-1	45%	





MDL = Method Detection Limit



RL = Reporting Limit

E = Indicates value exceeds calibration range

J = Indicates an estimated value

B = Indicates analyte found in associated method blank

N = Indicates presumptive evidence of a compound

SGS Accutest

Report of Analysis

Client Sample ID:	SB103-GWS
Lab Sample ID:	JC21103-8

AQ - Ground Water Matrix: Method:

SW846 8270D SW846 3510C Project:

05/26/16 Date Sampled: Date Received: 05/27/16 Percent Solids: n/a

Q

BMSMC, Building 5 Area, PR

Run #1	File ID Z111122.D	DF 1	Analyzed 05/28/16	By	Prep Date 05/27/16	Prep Batch OP94327	Analytical Batch EZ5550
Run #2							

Run #1 Run #2	Initial Volume 900 ml	Final Volume 1.0 ml			
17011 #2					

ABN TCL Special List

CAS No.	Compound	Result	RL	MDL	Units
95-57-8	2-Chlorophenol	ND	5.6	0.91	ug/l
59-50-7	4-Chloro-3-methyl phenol	ND	5.6	0.99	ug/l
120-83-2	2,4-Dichlorophenol	ND	2.2	1.4	ug/l
105-67-9	2,4-Dimethylphenol	ND	5.6	2.7	ug/l
51-28-5	2,4-Dinitrophenol	ND	11	1.7	ug/l
534-52-1	4,6-Dinitro-o-cresol	ND	5.6	1.4	ug/l
95-48-7	2-Methylphenol	ND	2.2	0.99	ug/l
	3&4-Methylphenol	ND	2.2	0.98	ug/l
88-75-5	2-Nitrophenol	ND	5.6	1.1	ug/l
100-02-7	4-Nitrophenol	ND	11	1.3	ug/l
87-86-5	Pentachlorophenol	ND	5.6	1.5	ug/l
108-95-2	Phenol	ND	2.2	0.44	ug/l
58-90-2	2,3,4,6-Tetrachlorophenol	ND	5.6	1.6	ug/l
95-95-4	2,4,5-Trichlorophenol	ND	5.6	1.5	ug/l
88-06-2	2,4,6-Trichlorophenol	ND	5.6	1.0	ug/l
83-32-9	Acenaphthene	ND	1.1	0.21	ug/l
208-96-8	Acenaphthylene	ND	1.1	0.15	ug/l
98-86-2	Acetophenone	ND	2.2	0.23	ug/l
120-12-7	Anthracene	ND	1.1	0.23	ug/l
1912-24-9	Atrazine	ND	2.2	0.50	ug/l
100-52-7	Benzaldehyde	ND	5.6	0.32	ug/l
56-55-3	Benzo(a)anthracene	ND	1.1	0.23	ug/l
50-32-8	Benzo(a)pyrene	ND	1.1	0.24	ug/l
205-99-2	Benzo(b)fluoranthene	ND	1.1	0.23	ug/l
191-24-2	Benzo(g,h,i)perylene	ND	1.1	0.38	ug/l
207-08-9	Benzo(k)fluoranthene	ND	1.1	0.23	ug/l
101-55-3	4-Bromophenyl phenyl ether	ND	2.2	0.45	ug/l
85-68-7	Butyl benzyl phthalate	ND	2.2	0.51	ug/l
92-52-4	1,1'-Biphenyl	ND	1.1	0.24	ug/l
91-58-7	2-Chloronaphthalene	ND	2.2	0.26	ug/l
106-47-8	4-Chloroaniline	ND	5.6	0.38	ug/l
86-74-8	Carhazole	ND	1.1	0.25	ug/l



ND = Not detected

MDL = Method Detection Limit

RL = Reporting Limit

E = Indicates value exceeds calibration range

J = Indicates an estimated value

B = Indicates analyte found in associated method blank

Client Sample ID: Lab Sample ID:

SB103-GWS JC21103-8

Matrix:

AQ - Ground Water

Method: Project:

SW846 8270D SW846 3510C BMSMC, Building 5 Area, PR

Date Sampled: Date Received:

05/26/16 05/27/16

Percent Solids: n/a

ABN TCL Special List

	•				
CAS No.	Compound	Result	RL	MDL	Units
105-60-2	Caprolactam	ND	2.2	0.72	ug/l
218-01-9	Chrysene	ND	1.1	0.20	ug/l
111-91-1	bis(2-Chloroethoxy)methane	ND	2.2	0.31	ug/l
111-44-4	bis(2-Chloroethyl)ether	ND	2.2	0.28	ug/l
108-60-1	bis(2-Chloroisopropyl)ether	ND	2.2	0.45	ug/l
7005-72-3	4-Chlorophenyl phenyl ether	ND	2.2	0.41	ug/l
121-14-2	2,4-Dinitrotoluene	ND	1.1	0.61	ug/l
606-20-2	2,6-Dinitrotoluene	ND	1.1	0.53	ug/l
91-94-1	3,3'-Dichlorobenzidine	ND	2.2	0.56	ug/l
53-70-3	Dibenzo(a,h)anthracene	ND	1.1	0.37	ug/l
132-64-9	Dibenzofuran	ND	5.6	0.24	ug/l
84-74-2	Di-n-butyl phthalate	ND	2.2	0.55	ug/l
117-84-0	Di-n-octyl phthalate	ND	2.2	0.26	ug/l
84-66-2	Diethyl phthalate	ND	2.2	0.29	ug/l
131-11-3	Dimethyl phthalate	ND	2.2	0.24	ug/l
117-81-7	bis(2-Ethylhexyl)phthalate	ND	2.2	1.8	ug/l
206-44-0	Fluoranthene	1.8	1.1	0.19	ug/l
86-73-7	Fluorene	ND	1.1	0.19	ug/l
118-74-1	Hexachlorobenzene	ND	1.1	0.36	ug/l
87-68-3	Hexachlorobutadiene	ND	1.1	0.55	ug/l
77-47-4	Hexachlorocyclopentadiene	ND	11	3.1	ug/l
67-72-1	Hexachloroethane	ND	2.2	0.43	ug/l
193-39-5	Indeno(1,2,3-cd)pyrene	ND	1.1	0.37	ug/l
78-59-1	Isophorone	ND	2.2	0.31	ug/l
90-12-0	1-Methylnaphthalene	ND	1.1	0.29	ug/l
91-57-6	2-Methylnaphthalene	ND	1.1	0.23	ug/l
88-74-4	2-Nitroaniline	ND	5.6	0.31	ug/l
99-09-2	3-Nitroaniline	ND	5.6	0.43	ug/l
100-01-6	4-Nitroaniline	ND	5.6	0.49	ug/l
98-95-3	Nitrobenzene	ND	2.2	0.71	ug/l
621-64-7	N-Nitroso-di-n-propylamine	ND	2.2	0.53	ug/l
86-30-6	N-Nitrosodiphenylamine	ND	5.6	0.25	ug/l
85-01-8	Phenanthrene	ND	1.1	0.19	ug/l
129-00-0	Pyrene	1.3	1.1	0.24	ug/l
95-94-3	1,2,4,5-Tetrachlorobenzene	ND	2.2	0.41	ug/l
CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Lim	its
367-12-4	2-Fluorophenol	50%		14-8	8%
4105.00.0	THE 1 15	-504			



ND = Not detected

4165-62-2

MDL = Method Detection Limit

37%

RL = Reporting Limit

E = Indicates value exceeds calibration range

Phenol-d5

J = Indicates an estimated value

10-110%

B = Indicates analyte found in associated method blank





Report of Analysis

Page 3 of 3

Client Sample ID: SB103-GWS

Lab Sample ID:

JC21103-8

Matrix: Method: Project:

AQ - Ground Water

SW846 8270D SW846 3510C BMSMC, Building 5 Area, PR

Date Sampled: 05/26/16 Date Received: 05/27/16

Percent Solids: n/a

ABN TCL Special List

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
118-79-6 4165-60-0	2,4,6-Tribromophenol Nitrobenzene-d5	74% 85%		39-149% 32-128%
321-60-8	2-Fluorobiphenyl	73%		35-119%
1718-51-0	Terphenyl-d14	58%		10-126%





MDL = Method Detection Limit

RL = Reporting Limit

J = Indicates an estimated value

B = Indicates analyte found in associated method blank



E = Indicates value exceeds calibration range

SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID: Lab Sample ID:

SB103-GWS

JC21103-8 AO - Ground Water Date Sampled: Date Received: 05/27/16

05/26/16

Matrix: Method:

SW846 8270D BY SIM SW846 3510C

Analyzed

06/01/16

Percent Solids: n/a

Project:

BMSMC, Building 5 Area, PR

DF

1

Prep Batch

Analytical Batch

Run #2

Run #1

LK 05/27/16

By

OP94327A

0

E4P880

Initial Volume Final Volume Run #1 900 ml 1.0 ml

4P16541.D

File ID

Run #2

CAS No.

91-20-3

123-91-1

CAS No.

4165-60-0

321-60-8

1718-51-0

Compound Naphthalene

1.4-Dioxane

Surrogate Recoveries

Nitrobenzene-d5

2-Fluorobiphenyl

Terphenyl-d14

Result ND

MDL 0.033

Prep Date

Units

ug/l ug/l

ND 0.11 0.054 Limits

Run#1 Run# 2

RL

0.11

82% 45% 60% 24-125% 19-127%

10-119%

tacl Infante Méndez IC = 1888

ND = Not detected RL = Reporting Limit MDL = Method Detection Limit

J = Indicates an estimated value B = Indicates analyte found in associated method blank



SGS Accutest

Report of Analysis

Page 1 of 1

Client Sample ID: SB103-GWS

Lab Sample ID:

JC21103-8

Matrix:

Method:

Project:

AQ - Ground Water

SW846-8015C (DAI)

BMSMC, Building 5 Area, PR

Date Sampled:

05/26/16 Date Received: 05/27/16

Percent Solids: n/a

Prep Date File ID DF Analyzed Prep Batch Analytical Batch By Run #1 GH105258.D 05/27/16 XPL GGH5303 1 n/a n/a Run #2

Low Molecular Alcohol List

Compound	Result	RL	MDL	Units	Q
Ethanol	ND	100	55	ug/l	
Isobutyl Alcohol	ND	100	36		
Isopropyl Alcohol	ND	100	68		
	ND	100	43		
	ND	100	87		
	ND	100	66	_	
Methanol	ND	200	71	ug/l	
Surrogate Recoveries	Run#1	Run# 2	Lim	its	
Hexanol	82%		56-1	45%	. e
Hexanol	87%		56-1	45%	SPE AS
	Ethanol Isobutyl Alcohol Isopropyl Alcohol n-Propyl Alcohol n-Butyl Alcohol sec-Butyl Alcohol Methanol Surrogate Recoveries Hexanol	Ethanol ND Isobutyl Alcohol ND Isopropyl Alcohol ND n-Propyl Alcohol ND n-Butyl Alcohol ND sec-Butyl Alcohol ND Methanol ND Surrogate Recoveries Run# 1 Hexanol 82%	Ethanol ND 100 Isobutyl Alcohol ND 100 Isopropyl Alcohol ND 100 n-Propyl Alcohol ND 100 n-Butyl Alcohol ND 100 sec-Butyl Alcohol ND 100 Methanol ND 200 Surrogate Recoveries Run#1 Run#2 Hexanol 82%	Ethanol	Ethanol





MDL = Method Detection Limit

RL = Reporting Limit

J = Indicates an estimated value

B = Indicates analyte found in associated method blank





E = Indicates value exceeds calibration range

Es	- Ca		and the	355		F(0.05)	AL.	FEST			TEL. 732-J		FAX:		3499/3	110		2730-	and the last	Marie a	Total Control	200	1000	TPE DE	100	Acceptance	1747	آد	21	Matrix Code
	7	. 14	الم	المعا) de	sa la	5	Project Name	110	Rel	***	40	Cerci	Oda a	erses	0		SCHOOL SECTION	NAME OF			81	2	1		T	No. of Concession, Name of Street, or other party of the Concession, Name of Street, or other pa	T	T	DW - Dreating Wa
Ė	ier.	1.1	- 6	I de	ν_{a}	~	-	Street		1CE	EUSL		7) 7 T C	// ·	MGT!	tac:	Section 1	200	4		180				1			1	GW - Ground Wa WW - Water
		Wę		16	5161	Zφ	-	Cay	_		AM	Conpa	andormani Ny Indrina	-11-	to-ont (rum R	eport be	4	_	8270	015	00	Es		-)					SW - Surface Wa SQ - Soil SL Skidge
U	C	as	2	N.	Υ_	E-mail		HU	MO	Ca o	PR	Street A				33						3					1	1		SED-Sedemen
Ţ	201	y ·	Ta	410	25			- Industry												Men	0	1	102				1			LIQ - Other Liqui ARI - Air SQL - Other Set
Ď	14	1-7	15	ť-,	041	00		Clerel Punchs	obse Oro	w 1		Cay				4386		24		3	3		200			1			8	VrP - Wipe FB-Feld Blank
R	W.	,7	Tay	ler,	D.L	ndst	7347	TOPIC SALES	actor*			Attento	•		77					1	Z	8	Q							EB-Equipment (b) HB-Throne (b) TB-Tro-(b)ann
							T		F	Q.	Calterton				H	T		i.T	To I	100		4	Þξ							
3	F	eld ID/	Point o	d Com	ection	115	-	MCCHOI VAL		-45	1-	57	Man.	-	9	ğ ğ	100 A	10 Physical March	DOM:	1	1	20	53	. 1	-1	1	1			LAB USE ONL
1	-	47	G	-	- 7 -		4	15/20	1	05/26/	1535	TY	6W	6	3		3		П	X	7	X	X		-					K10
	E	05	25	2	016	-	1	_	0	5125/K	1548	NK	EB	6	3	Н	9		П	X	1	X	X							V238
	5	5/6	23	(Ζ,	- 5	3 4	24	-	10	5/24/14	1030	11	50	2	₩	Н	2	Н	11	1X	X	_	X	-	-	\perp	\perp	1	N.	Dro
+	킷	310	3	ų	2	7 5	2/	_	0	2/14/16	1030	11	50	2	₩	H	2	1	Н	1X	7-		Ķ	-	1	+	+	+-	W.	H3T+
	À.	7/	4	ι	_	11	-		- 12	724//5	1035	11	50	2	++	1	2	1	H	1X	7		X	+	-	+	1	-	\vdash	F53
	17	0.17	3	-4	=W	K	-		24	3 /26/16 Clay / 1	1152	M	50	14		H	2	H	+	1	4	_	7	+	-	+	-	-		
+	-	RIC	3	_	- 1. 1	2	_	-	0	<u> </u>	1230	-	GW	2	13	Н	2	Η-	Н	X	^	_	_	-	1	+	+-	-		-
-	#i	310	3-	9	LHV	2	+		1	इंग्याह	1400	TT	ΘW.	5	13	Н	2	H	Н	X	X.	-	K	+	-	+	+	-	5	Ac
1	-		5				\neg	-		-	-			_	tt	H	Н	-	H	1			-	+	HITH	4-4	SESSI	AENT	14	76
				100			1			-					Ħ	\forall	П		\vdash					+	LABE	VE	FIFIC	TIO	_4	=
	Š.				_			-							П	П						12.1								-
9	_	Turnamo	_						03 A	A DAZ	7.34	-	Commerc	_			ratorma		P Cate	B-7+	1160	X	A A	A C	OCCUPANT OF	/ Spe	cini limate	chore	ESS.	CT STORY TON
3	Sant.	10 Beek	ota De	rfo	r 5o	يهكر إر	rple	<u> </u>		_			Carrange			-		MYAS	P Categ			n	10 10	Kep	brT:	1-	Ach	11/4	2 P/V	Halene D
i] 35	ry Rush					-			= 1		 	PULLT† (NJ Rođes	Limpi 3+ ad	4)		H	Store				W	th	SV	OC !	Me	thod	2:	270	D
Į		y RUSH w Strike	Gal	- A4	7060	1) Sau 145 Sa	240	25					Commerc	mage.				Other					0-1-	17/2000				-		
1	j						- 1		-	=		Comme	74.7 (AT 1800)	ne Alaban Republis Di			2000 (Pig gi 1911 (c)	pering Results	-005	uristrary	Ť	-	_	-	_	_	_	_	_	
ź	POR S	Rush To	444	1000	VA Lab	int.	-	-	Samol	Codes	and from all con-	NJ Page	cod = Res	met + QC	Summ	ey - P	what Po	in dete		_			de inve	nitory is	venñe	d upo	n reces	st in In	e Labo	restory
			_	(Congress)	No. of Contract of			170									ую рос	4404.10	M, INCH	wane c	ourier e	September	w.			ALC: NO.				

JC21103: Chain of Custody Page 1 of 2

EXECUTIVE NARRATIVE

SDG No:

JC21103

Laboratory:

Accutest, New Jersey

Analysis:

SW846-8270D

Number of Samples:

8

Location:

BMSMC, Building 5 Area

Humacao, PR

SUMMARY: Eight (8) samples were analyzed for the ABN TCL list following method SW846-8270D; Naphthalene and 1,4-Dioxane were also analyzed by SW846-8270D using the selective ion monitoring (SIM) technique. The sample results were assessed according to USEPA data validation guidance documents in the following order of precedence: EPA Hazardous Waste Support Section, SOP HW-35A, July 2015 –Revision 0. Semivolatile Data Validation. The QC criteria and data validation actions listed on the data review worksheets are from the primary guidance document, unless otherwise noted.

Results are valid and can be used for decision making purposes.

Critical issues:

None

Major:

None

Minor:

None

Critical findings: Major findings:

None None

Minor findings:

1. Initial and continuing calibration verifications meet the required criteria. Analytes not meeting the method % difference criteria meet the guidance document performance criteria for continuing calibration verification of +25 or 40 %, no action taken.

Analytes not meeting the % difference continuing calibration criteria were qualified as estimated (J) or (UJ) in affected samples.

No closing calibration verification included in data package. No action taken, professional judgment.

2. MS/MSD % recoveries and RPD outside the in-house limits for several analytes in the MS/MSD QC sample for this sample batch. No action taken, MS/MSD results apply to unspiked sample. Unspiked samples were from another project.

No aqueous matrix MS/MSD sample was analyzed with this sample batch. Blank spike/blank spike duplicate used to assess accuracy. MS/MSD % recoveries and PRP within laboratory control limits. No action taken, professional judgment.

COMMENTS:

Results are valid and can be used for decision making purposes.

Reviewers Name:

Rafael Infante

Chemist License 1888

Signature:

Date:

Wafaul arfaut

SAMPLE ORGANIC DATA SAMPLE SUMMARY

Sample ID: JC21103-1

Sample location: BMSMC Building 5 Area

Sampling date: 5/25/2016 Matrix: Groundwater

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
2-Chlorophenol	5.4	ug/l	1	-	U	Yes
4-Chloro-3-methyl phenol	5.4	ug/l	1	-	U	Yes
2,4-Dichlorophenol	2.2	ug/i	1	-	U	Yes
2,4-Dimethylphenol	5.4	ug/l	1	-	U	Yes
2,4-Dinitrophenol	11	ug/l	1	-	U	Yes
4,6-Dinitro-o-cresol	5.4	ug/l	1	-	U	Yes
2-Methylphenol	2.2	ug/l	1	-	U	Yes
3&4-Methylphenol	2.2	ug/l	1	-	U	Yes
2-Nitrophenol	5.4	ug/l	1	-	U	Yes
4-Nitrophenol	11	ug/l	1	-	U	Yes
Pentachlorophenol	5.4	ug/l	1	-	U	Yes
Phenol	2.2	ug/l	1	-	UJ	Yes
2,3,4,6-Tetrachlorophenol	5.4	ug/l	1	-	U	Yes
2,4,5-Trichlorophenol	5.4	ug/l	1	-	U	Yes
2,4,6-Trichlorophenol	5.4	ug/i	1	-	U	Yes
Acenaphthene	1.1	ug/l	1	-	U	Yes
Acenaphthylene	1.1	ug/l	1	-	U	Yes
Acetophenone	2.2	ug/l	1	-	U	Yes
Anthracene	1.1	ug/l	1	-	U	Yes
Atrazine	2.2	ug/l	1	-	U	Yes
Benzaldehyde	5.4	ug/l	1	-	U	Yes
Benzo(a) anthracene	1.1	ug/l	1	-	U	Yes
Benzo(a)pyrene	1.1	ug/l	1	-	U	Yes
Benzo(b)fluoranthene	1.1	ug/l	1	-	U	Yes
Benzo(g,h,i)perylene	1.1	ug/l	1	-	U	Yes
Benzo(k)fluoranthene	1.1	ug/l	1	-	U	Yes
4-Bromophenyl phenyl ether	1.1	ug/l	1	-	U	Yes
Butyl benzyl phthalate	2.2	ug/l	1	-	U	Yes
1,1'-Biphenyl	1.1	ug/l	1	-	U	Yes
2-Chloronaphthalene	2.2	ug/l	1	-	U	Yes
4-Chloroaniline	5.4	ug/l	1	-	U	Yes
Carbazole	1.1	ug/l	1	-	U	Yes
Caprolactam	2.2	ug/l	1	-	U	Yes
Chrysene	1.1	ug/l	1	-	U	Yes
bis(2-Chloroethoxy)methane	2.2	ug/l	1	-	U	Yes
bis (2-Chloroethyl) ether	2.2	ug/l	1	-	UJ	Yes

bis(2-Chloroisopropyl)ether	2.2	ug/l	1	528	U	Yes
4-Chlorophenyl phenyl ether	2.2	ug/l	1	(10)	U	Yes
2,4-Dinitrotoluene	1.1	ug/l	1	V. T.	U	Yes
2,6-Dinitrotoluene	1.1	ug/l	1	-	U	Yes
3,3'-Dichlorobenzidine	2.2	ug/l	1		U	Yes
1,4-Dioxane	6.6	ug/l	1	-	_	Yes
Dibenzo(a,h)anthracene	1.1	ug/l	1	-	U	Yes
Dibenzofuran	5.4	ug/l	1	-	U	Yes
Di-n-butyl phthalate	2.2	ug/l	1	121	U	Yes
Di-n-octyl phthalate	2.2	ug/l	1	-	U	Yes
Diethyl phthalate	2.2	ug/l	1		U	Yes
Dimethyl phthalate	2.2	ug/l	1	-	U	Yes
bis (2-Ethylhexyl) phthalate	2.2	ug/l	1	-	U	Yes
Fluoranthene	1.1	ug/l	1	-	U	Yes
Fluorene	1.1	ug/l	1	140	U	Yes
Hexachlorobenzene	1.1	ug/l	1	-	U	Yes
Hexachlorobutadiene	1.1	ug/l	1	-	U	Yes
Hexachlorocyclopentadiene	11	ug/l	1	-	U	Yes
Hexachloroethane	2.2	ug/l	1	-	U	Yes
Indeno(1,2,3-cd)pyrene	1.1	ug/l	1	-	U	Yes
Isophorone	2.2	ug/l	1	-	Ų	Yes
1-Methylnaphthalene	1.1	ug/l	1	-	U	Yes
2-Methylnaphthalene	1.1	ug/l	1	-	U	Yes
2-Nitroaniline	5.4	ug/l	1	-	UJ	Yes
3-Nitroaniline	5.4	ug/l	1	_	UJ	Yes
4-Nitroaniline	5.4	ug/l	1	-	U	Yes
Nitrobenzene	2.2	ug/l	1	7	U	Yes
N-Nitroso-di-n-propylamine	2.2	ug/l	1	-	LU.	Yes
Nitrosodiphenylamine	5.4	ug/l	1	-	Ų	Yes
Phenanthrene	1.1	ug/i	1	-	U	Yes
Pyrene	1.1	ug/l	1	-	U	Yes
1,2,4,5-Tetrachlorobenzene	2.2	ug/l	1	-	U	Yes
METHOD:	8270D (SI	M)				
Naphthalene	0.11	ug/l	1	(5.)	U	Yes

Sample location: BMSMC Building 5 Area

Sampling date: 5/25/2016

Matrix: AQ - Equipment Blank

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
2-Chlorophenol	5.0	ug/l	1	-	U	Yes
4-Chloro-3-methyl phenol	5.0	ug/l	1	-	U	Yes
2,4-Dichlorophenol	2.0	ug/l	1	-	U	Yes
2,4-Dimethylphenol	5.0	ug/l	1	-	U	Yes
2,4-Dinitrophenol	10	ug/l	1	-	U	Yes
4,6-Dinitro-o-cresol	5.0	ug/l	1	-	U	Yes
2-Methylphenol	2.0	ug/l	1	-	U	Yes
3&4-Methylphenol	2.0	ug/l	1	-	U	Yes
2-Nitrophenol	5.0	ug/l	1	-	U	Yes
4-Nitrophenol	10	ug/l	1	-	U	Yes
Pentachlorophenol	5.0	ug/l	1	-	U	Yes
Phenol	2.0	ug/l	1	-	U	Yes
2,3,4,6-Tetrachlorophenol	5.0	ug/l	1	-	U	Yes
2,4,5-Trichlorophenol	5.0	ug/l	1	-	U	Yes
2,4,6-Trichlorophenol	5.0	ug/l	1	-	U	Yes
Acenaphthene	1.0	ug/l	1	-	U	Yes
Acenaphthylene	1.0	ug/l	1	-	U	Yes
Acetophenone	2.0	ug/i	1	-	U	Yes
Anthracene	1.0	ug/l	1	-	U	Yes
Atrazine	2.0	ug/l	1	-	U	Yes
Benzaldehyde	5.0	ug/l	1	-	U	Yes
Benzo(a)anthracene	1.0	ug/l	1	-	U	Yes
Benzo(a)pyrene	1.0	ug/l	1	-	U	Yes
Benzo(b)fluoranthene	1.0	ug/l	1	•	U	Yes
Benzo(g,h,i)perylene	1.0	ug/l	1	-	U	Yes
Benzo(k)fluoranthene	1.0	ug/l	1	-	U	Yes
4-Bromophenyl phenyl ether	1.0	ug/l	1	-	U	Yes
Butyl benzyl phthalate	2.0	ug/l	1	-	U	Yes
1,1'-Biphenyl	1.0	ug/l	1	-	Ų	Yes
2-Chloronaphthalene	2.0	ug/l	1	•	U	Yes
4-Chloroaniline	5.0	ug/l	1	-	U	Yes
Carbazole	1.0	ug/l	1	-	U	Yes
Caprolactam	2.0	ug/l	1	~	U	Yes
Chrysene	1.0	ug/l	1	-	U	Yes
bis(2-Chloroethoxy)methane	2.0	ug/l	1	-	U	Yes
bis (2-Chloroethyl) ether	2.0	ug/l	1	-	UJ	Yes
bis (2-Chlorois opropyl) ether	2.0	ug/l	1	-	U	Yes

4-Chlorophenyl phenyl ether	2.0	ug/l	1	-	U	Yes
2,4-Dinitrotoluene	1.0	ug/l	1	-	U	Yes
2,6-Dinitrotoluene	1.0	ug/l	1	-	U	Yes
3,3'-Dichlorobenzidine	2.0	ug/l	1	12	U	Yes
Dibenzo(a,h)anthracene	1.0	ug/l	1	-	U	Yes
Dibenzofuran	5.0	ug/l	1		U	Yes
Di-n-butyl phthalate	2.0	ug/l	1	-	U	Yes
Di-n-octyl phthalate	2.0	ug/l	1		U	Yes
Diethyl phthalate	2.0	ug/l	1	-	U	Yes
Dimethyl phthalate	2.0	ug/l	1	-	U	Yes
bis(2-Ethylhexyl)phthalate	2.0	ug/l	1	-	U	Yes
Fluoranthene	1.0	ug/l	1	_	U	Yes
Fluorene	1.0	ug/l	1	-	U	Yes
Hexachlorobenzene	1.0	ug/l	1	-	U	Yes
Hexachlorobutadiene	1.0	ug/l	1	-	U	Yes
Hexachlorocyclopentadiene	10	ug/l	1	-	U	Yes
Hexachloroethane	2.0	ug/l	1	-	U	Yes
Indeno(1,2,3-cd)pyrene	1.0	ug/l	1		U	Yes
Isophorone	2.0	ug/l	1		υ	Yes
1-Methylnaphthalene	1.0	ug/l	1	-	U	Yes
2-Methylnaphthalene	1.0	ug/l	1	-	U	Yes
2-Nitroaniline	5.0	ug/l	1	-	UJ	Yes
3-Nitroaniline	5.0	ug/l	1	-	UJ	Yes
4-Nitroaniline	5.0	ug/l	1	-	U	Yes
Nitrobenzene	2.0	ug/l	1	-	U	Yes
N-Nitroso-di-n-propylamine	2.0	ug/l	1	-	UJ	Yes
Nitrosodiphenylamine	5.0	ug/l	1		U	Yes
Phenanthrene	1.0	ug/l	1	-	U	Yes
Pyrene	1.0	ug/l	1	-	U	Yes
1,2,4,5-Tetrachlorobenzene	2.0	ug/l	1	-	U	Yes
METHOD:	8270D (SI	M)				
Naphthalene	0.10	ug/l	1		U	Yes
1,4-Dioxane	0.10	ug/l	1		U	Yes
A) T WICHUIG	0.10	u6/ i	1	4773	U	1 5

Sample location: BMSMC Building 5 Area

Sampling date: 5/26/2016

Matrix: Soil

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
2-Chlorophenol	73	ug/kg	1	-	U	Yes
4-Chloro-3-methyl phenol	180	ug/kg	1	-	U	Yes
2,4-Dichlorophenol	180	ug/kg	1	-	U	Yes
2,4-Dimethylphenol	180	ug/kg	1	-	U	Yes
2,4-Dinitrophenol	180	ug/kg	1	-	U	Yes
4,6-Dinitro-o-cresol	180	ug/kg	1	-	U	Yes
2-Methylphenol	73	ug/kg	1	-	U	Yes
3&4-Methylphenol	73	ug/kg	1	-	U	Yes
2-Nitrophenol	180	ug/kg	1	-	U	Yes
4-Nitrophenol	370	ug/kg	1	-	UJ	Yes
Pentachlorophenol	180	ug/kg	1	-	UJ	Yes
Phenol	73	ug/kg	1	-	U	Yes
2,3,4,6-Tetrachlorophenol	180	ug/kg	1	-	U	Yes
2,4,5-Trichlorophenol	180	ug/kg	1	-	U	Yes
2,4,6-Trichlorophenol	180	ug/kg	1	-	U	Yes
Acenaphthene	37	ug/kg	1	-	U	Yes
Acenaphthylene	37	ug/kg	1	-	U	Yes
Acetophenone	180	ug/kg	1	-	U	Yes
Anthracene	37	ug/kg	1	-	U	Yes
Atrazine	73	ug/kg	1	-	U	Yes
Benzo(a)anthracene	37	ug/kg	1	-	U	Yes
Benzo(a)pyrene	37	ug/kg	1	-	U	Yes
Benzo(b)fluoranthene	37	ug/kg	1	-	U	Yes
Benzo(g,h,i)perylene	37	ug/kg	1	-	U	Yes
Benzo(k)fluoranthene	37	ug/kg	1	-	U	Yes
4-Bromophenyl phenyl ether	73	ug/kg	1	-	U	Yes
Butyl benzyl phthalate	73	ug/kg	1	-	U	Yes
1,1'-Biphenyl	73	ug/kg	1	-	U	Yes
Benzaldehyde	180	ug/kg	1	-	U	Yes
2-Chloronaphthalene	73	ug/kg	1	-	U	Yes
4-Chloroaniline	180	ug/kg	1	-	U	Yes
Carbazole	73	ug/kg	1	-	U	Yes
Caprolactam	73	ug/kg	1	-	U	Yes
Chrysene	37	ug/kg	1	-	U	Yes
bis(2-Chloroethoxy)methane	73	ug/kg	1	-	U	Yes
bis(2-Chloroethyl)ether	73	ug/kg	1	-	U	Yes
bis(2-Chloroisopropyl)ether	73	ug/kg	1	-	U	Yes

4-Chlorophenyl phenyl ether	73	ug/kg	1	82	U	Yes
2,4-Dinitrotoluene	37	ug/kg	1	-	U	Yes
2,6-Dinitrotoluene	37	ug/kg	1	7	U	Yes
3,3'-Dichlorobenzidine	73	ug/kg	1	-	U	Yes
Dibenzo(a,h)anthracene	37	ug/kg	1	-	U	Yes
Dibenzofuran	73	ug/kg	1	-	U	Yes
Di-n-butyl phthalate	73	ug/kg	1	-	U	Yes
Di-n-octyl phthalate	73	ug/kg	1	7.	U	Yes
Diethyl phthalate	73	ug/kg	1	2	U	Yes
Dimethyl phthalate	73	ug/kg	1	-	U	Yes
bis (2-Ethylhexyl) phthalate	73	ug/kg	1	1.7	U	Yes
Fluoranthene	37	ug/kg	1	-	U	Yes
Fluorene	37	ug/kg	1	-	U	Yes
Hexachlorobenzene	73	ug/kg	1	-	U	Yes
Hexachlorobutadiene	37	ug/kg	1	-	U	Yes
Hexachlorocyclopentadiene	410	ug/kg	1	-	U	Yes
Hexachloroethane	180	ug/kg	1	- 2	U	Yes
Indeno(1,2,3-cd)pyrene	37	ug/kg	1	-	U	Yes
Isophorone	73	ug/kg	1		U	Yes
1-Methylnaphthalene	73	ug/kg	1	-	U	Yes
2-Methylnaphthalene	73	ug/kg	1	(5)	U	Yes
2-Nitroaniline	180	ug/kg	1	-	U	Yes
3-Nitroaniline	180	ug/kg	1	-	U	Yes
4-Nitroaniline	180	ug/kg	1	-	U	Yes
Nitrobenzene	73	ug/kg	1	-	U	Yes
N-Nitroso-di-n-propylamine	73	ug/kg	1	-	U	Yes
Nitrosodiphenylamine	180	ug/kg	1	-	U	Yes
Phenanthrene	37	ug/kg	1	-	U	Yes
Pyrene	37	ug/kg	1	-	U	Yes
1,2,4,5-Tetrachlorobenzene	180	ug/kg	1	-	U	Yes

METHOD: 8270D (SIM)

Naphthalene	3.7	ug/kg	1	-	U	Yes
1,4-Dioxane	3.7	ug/kg	1		U	Yes

Sample location: BMSMC Building 5 Area

Sampling date: 5/26/2016

Matrix: Soil

METHOD:	82/00					
Analyte Name	Result		Dilution Factor	Lab Flag	Validation	Reportable
2-Chlorophenol	70	ug/kg	1	-	U	Yes
4-Chloro-3-methyl phenol	170	ug/kg	1	-	U	Yes
2,4-Dichlorophenol	170	ug/kg	1	-	U	Yes
2,4-Dimethylphenol	170	ug/kg	1	-	U	Yes
2,4-Dinitrophenol	170	ug/kg	1	-	U	Yes
4,6-Dinitro-o-cresol	170	ug/kg	1	-	U	Yes
2-Methylphenol	70	ug/kg	1	-	U	Yes
3&4-Methylphenol	70	ug/kg	1	-	U	Yes
2-Nitrophenol	170	ug/kg	1	-	U	Yes
4-Nitrophenol	350	ug/kg	1	-	U	Yes
Pentachlorophenol	170	ug/kg	1	-	U	Yes
Phenol	70	ug/kg	1	-	U	Yes
2,3,4,6-Tetrachlorophenol	170	ug/kg	1	-	U	Yes
2,4,5-Trichlorophenol	170	ug/kg	1 🧓	•	U	Yes
2,4,6-Trichlorophenol	170	ug/kg	1	-	U	Yes
Acenaphthene	35	ug/kg	1	-	U	Yes
Acenaphthylene	35	ug/kg	1	-	U	Yes
Acetophenone	170	ug/kg	1	-	U	Yes
Anthracene	35	ug/kg	1	-	U	Yes
Atrazine	70	ug/kg	1	-	U	Yes
Benzo(a)anthracene	35	ug/kg	1	-	U	Yes
Benzo(a)pyrene	35	ug/kg	1	-	U	Yes
Benzo(b)fluoranthene	35	ug/kg	1	-	U	Yes
Benzo(g,h,i)perylene	35	ug/kg	1	-	U	Yes
Benzo(k)fluoranthene	35	ug/kg	1	-	U	Yes
4-Bromophenyl phenyl ether	70	ug/kg	1		U	Yes
Butyl benzyi phthalate	70	ug/kg	1	-	UJ	Yes
1,1'-Biphenyl	70	ug/kg	1	-	U	Yes
Benzaldehyde	170	ug/kg	1	•	U	Yes
2-Chloronaphthalene	70	ug/kg	1	-	U	Yes
4-Chloroaniline	170	ug/kg	1	-	U	Yes
Carbazole	70	ug/kg	1	-	U	Yes
Caprolactam	70	ug/kg	1	-	U	Yes
Chrysene	35	ug/kg	1	-	U	Yes
bis(2-Chloroethoxy)methane	70	ug/kg	1	-	U	Yes
bis(2-Chloroethyl)ether	70	ug/kg	1	-	Ų	Yes
bis(2-Chloroisopropyl)ether	70	ug/kg	1	-	U	Yes

4-Chlorophenyl phen	yl ether 70	ug/kg	1	-	U	Yes
2,4-Dinitrotoluene	35	ug/kg	1	-	U	Yes
2,6-Dinitrotoluene	35	ug/kg	1	_	U	Yes
3,3'-Dichlorobenzidin	e 70	ug/kg	1	-	U	Yes
Dibenzo(a,h)anthrace	ene 35	ug/kg	1	-	U	Yes
Dibenzofuran	70	ug/kg	1	-	U	Yes
Di-n-butyl phthalate	70	ug/kg	1	-	U	Yes
Di-n-octyl phthalate	70	ug/kg	1	-	U	Yes
Diethyl phthalate	70	ug/kg	1	5-1	U	Yes
Dimethyl phthalate	70	ug/kg	1	-	U	Yes
bis (2-Ethylhexyl) phth	alate 70	ug/kg	1	-	U	Yes
Fluoranthene	35	ug/kg	1	-	U	Yes
Fluorene	35	ug/kg	1	-	U	Yes
Hexachlorobenzene	70	ug/kg	1	_	U	Yes
Hexachlorobutadiene	35	ug/kg	1	-	U	Yes
Hexachlorocyclopenta	adiene 350	ug/kg	1		U	Yes
Hexachloroethane	170	ug/kg	1	-	U	Yes
Indeno(1,2,3-cd)pyrei	ne 35	ug/kg	1	-	U	Yes
Isophorone	70	ug/kg	1	-	U	Yes
1-Methylnaphthalene	2 70	ug/kg	1	-	U	Yes
2-Methylnaphthalene	. 70	ug/kg	1	-	U	Yes
2-Nitroaniline	170	ug/kg	1	-	UJ	Yes
3-Nitroaniline	170	ug/kg	1	-	U	Yes
4-Nitroaniline	170	ug/kg	1	-	U	Yes
Nitrobenzene	70	ug/kg	1	-	U	Yes
N-Nitroso-di-n-propyl		ug/kg	1	-	U	Yes
Nitrosodiphenylamine	e 170	ug/kg	1	-	U	Yes
Phenanthrene	35	ug/kg	1	-	U	Yes
Pyrene	35	ug/kg	1		U	Yes
1,2,4,5-Tetrachlorobe	nzene 170	ug/kg	1	12	U	Yes
S						
	METHOD: 8270D (•				
Naphthalene	3.6	ug/kg	1	-	U	Yes
1,4-Dioxane	3.6	ug/kg	1	-	U	Yes

Sample location: BMSMC Building 5 Area

Sampling date: 5/26/2016

Matrix: Soil

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
2-Chlorophenol	77	ug/kg	1	-	U	Yes
4-Chloro-3-methyl phenol	190	ug/kg	1	-	U	Yes
2,4-Dichlorophenol	190	ug/kg	1	-	U	Yes
2,4-Dimethylphenol	190	ug/kg	1	-	U	Yes
2,4-Dinitrophenol	190	ug/kg	1	-	U	Yes
4,6-Dinitro-o-cresol	190	ug/kg	1	-	U	Yes
2-Methylphenol	77	ug/kg	1	-	U	Yes
3&4-Methylphenol	77	ug/kg	1	-	U	Yes
2-Nitrophenol	190	ug/kg	1	-	UJ	Yes
4-Nitrophenol	350	ug/kg	1	-	U	Yes
Pentachlorophenol	190	ug/kg	1	-	U	Yes
Phenol	77	ug/kg	1	-	U	Yes
2,3,4,6-Tetrachlorophenol	190	ug/kg	1	-	UJ	Yes
2,4,5-Trichlorophenol	190	ug/kg	1	-	U	Yes
2,4,6-Trichlorophenol	190	ug/kg	1	-	U	Yes
Acenaphthene	38	ug/kg	1	-	U	Yes
Acenaphthylene	38	ug/kg	1	-	U	Yes
Acetophenone	190	ug/kg	1	-	U	Yes
Anthracene	56.8	ug/kg	1	-	-	Yes
Atrazine	77	ug/kg	1	-	U	Yes
Benzo(a)anthracene	76.3	ug/kg	1	-	-	Yes
Benzo(a)pyrene	48.7	ug/kg	1	-	-	Yes
Benzo(b)fluoranthene	66.4	ug/kg	1	-	-	Yes
Benzo(g,h,i)perylene	19.8	ug/kg	1	J	UJ	Yes
Benzo(k)fluoranthene	24.9	ug/kg	1	J	UJ	Yes
4-Bromophenyl phenyl ether	77	ug/kg	1	-	U	Yes
Butyl benzyl phthalate	77	ug/kg	1	-	UJ	Yes
1,1'-Biphenyl	77	ug/kg	1	-	U	Yes
Benzaldehyde	190	ug/kg	1	-	U	Yes
2-Chloronaphthalene	7 7	ug/kg	1	-	U	Yes
4-Chloroaniline	190	ug/kg	1	-	U	Yes
Carbazole	77	ug/kg	1	-	U	Yes
Caprolactam	77	ug/kg	1	-	U	Yes
Chrysene	25.9	ug/kg	1	J	UJ	Yes
bis(2-Chloroethoxy)methane	77	ug/kg	1	-	U	Yes
bis(2-Chloroethyl)ether	77	ug/kg	1	-	U	Yes
bis(2-Chloroisopropyl)ether	77	ug/kg	1	-	U	Yes
4-Chlorophenyl phenyl ether	77	ug/kg	1	-	U	Yes

2,4-Dinitrotoluene	38	ug/kg	1		U	Yes
2,6-Dinitrotoluene	38	ug/kg	1	-	UJ	Yes
3,3'-Dichlorobenzidine	77	ug/kg	1	-	U	Yes
Dibenzo(a,h)anthracene	38	ug/kg	1		U	Yes
Dibenzofuran	77	ug/kg	1		U	Yes
Di-n-butyl phthalate	77	ug/kg	1		U	Yes
Di-n-octyl phthalate	77	ug/kg	1		U	Yes
Diethyl phthalate	77	ug/kg	1	-	U	Yes
Dimethyl phthalate	77	ug/kg	1	-	U	Yes
bis(2-Ethylhexyl)phthalate	77	ug/kg	1	3.40	U	Yes
Fluoranthene	501	ug/kg	1	-	-	Yes
Fluorene	38	ug/kg	1		U	Yes
Hexachlorobenzene	77	ug/kg	1	-	U	Yes
Hexachlorobutadiene	38	ug/kg	1	_	UJ	Yes
Hexachlorocyclopentadiene	380	ug/kg	1	(-)	U	Yes
Hexachloroethane	190	ug/kg	1	-	U	Yes
Indeno(1,2,3-cd)pyrene	21.2	ug/kg	1	1	UJ	Yes
Isophorone	77	ug/kg	1	-	U	Yes
1-Methylnaphthalene	77	ug/kg	1	-	U	Yes
2-Methylnaphthalene	77	ug/kg	1		U	Yes
2-Nitroaniline	190	ug/kg	1	7.0	נט	Yes
3-Nitroaniline	190	ug/kg	1	2	UJ	Yes
4-Nitroaniline	190	ug/kg	1	-	U	Yes
Nitrobenzene	77	ug/kg	1	7.5	U	Yes
N-Nitroso-di-n-propylamine	77	ug/kg	1	-	U	Yes
Nitrosodiphenylamine	190	ug/kg	1	-	U	Yes
Phenanthrene	38	ug/kg	1	-	U	Yes
Pyrene	321	ug/kg	1	-	-	Yes
1,2,4,5-Tetrachlorobenzene	190	ug/kg	1	•	U	Yes
METHOD:	8270D (S	IM)				
Naphthalene	3.8	ug/kg	1	-	U	Yes
1,4-Dioxane	3.8	ug/kg	1	-	U	Yes

Sample location: BMSMC Building 5 Area

Sampling date: 5/26/2016

Matrix: Soil

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
2-Chlorophenol	80	ug/kg	1	-	U	Yes
4-Chloro-3-methyl phenol	200	ug/kg	1	-	U	Yes
2,4-Dichlorophenol	200	ug/kg	1	-	U	Yes
2,4-Dimethylphenol	200	ug/kg	1	-	U	Yes
2,4-Dinitrophenol	200	ug/kg	1	-	U	Yes
4,6-Dinitro-o-cresol	200	ug/kg	1	-	U	Yes
2-Methylphenol	80	ug/kg	1	-	U	Yes
3&4-Methylphenol	80	ug/kg	1	-	U	Yes
2-Nitrophenol	200	ug/kg	1	-	U	Yes
4-Nitrophenol	400	ug/kg	1	-	U	Yes
Pentachlorophenol	200	ug/kg	1	-	U	Yes
Phenol	80	ug/kg	1	-	U	Yes
2,3,4,6-Tetrachlorophenol	200	ug/kg	1	-	U	Yes
2,4,5-Trichlorophenol	200	ug/kg	1	-	U	Yes
2,4,6-Trichlorophenol	200	ug/kg	1	-	U	Yes
Acenaphthene	40	ug/kg	1	-	U	Yes
Acenaphthylene	40	ug/kg	1	-	U	Yes
Acetophenone	200	ug/kg	1	-	U	Yes
Anthracene	40	ug/kg	1	-	U	Yes
Atrazine	80	ug/kg	1	-	U	Yes
Benzo(a)anthracene	40	ug/kg	1	-	U	Yes
Benzo(a)pyrene	40	ug/kg	1	-	U	Yes
Benzo(b)fluoranthene	40	ug/kg	1	-	U	Yes
Benzo(g,h,i)perylene	40	ug/kg	1	-	U	Yes
Benzo(k)fluoranthene	40	ug/kg	1	-	U	Yes
4-Bromophenyl phenyl ether	80	ug/kg	1	-	U	Yes
Butyl benzyl phthalate	80	ug/kg	1	-	UJ	Yes
1,1'-Biphenyl	80	ug/kg	1	-	U	Yes
Benzaldehyde	200	ug/kg	1	-	U	Yes
2-Chloronaphthalene	80	ug/kg	1	-	U	Yes
4-Chloroaniline	200	ug/kg	1	-	U	Yes
Carbazole	80	ug/kg	1	-	U	Yes
Caprolactam	80	ug/kg	1	-	U	Yes
Chrysene	40	ug/kg	1	-	U	Yes
bis(2-Chloroethoxy)methane	80	ug/kg	1	-	U	Yes
bis(2-Chloroethyl)ether	80	ug/kg	1	-	U	Yes
bis (2-Chloro is opropyl) ether	80	ug/kg	1	-	U	Yes

4-Chlorophenyl phenyl ether	80	ug/kg	1	2	U	Yes
2,4-Dinitrotoluene	40	ug/kg	1		Ü	Yes
2,6-Dinitrotoluene	40	ug/kg	1		Ü	Yes
3,3'-Dichlorobenzidine	80	ug/kg	1		Ü	Yes
Dibenzo(a,h)anthracene	40	ug/kg	1		Ü	Yes
Dibenzofuran	80	ug/kg	1		Ü	Yes
Di-n-butyl phthalate	80	ug/kg	1	*	Ü	Yes
Di-n-octyl phthalate	80	ug/kg	1		Ü	Yes
Diethyl phthalate	80	ug/kg	1	_	Ū	Yes
Dimethyl phthalate	80	ug/kg	1		Ū	Yes
bis(2-Ethylhexyl)phthalate	80	ug/kg	1		U	Yes
Fluoranthene	40	ug/kg	1		U	Yes
Fluorene	40	ug/kg	1		U	Yes
Hexachlorobenzene	80	ug/kg	1		U	Yes
Hexachlorobutadiene	40	ug/kg	1	+	U	Yes
Hexachlorocyclopentadiene	400	ug/kg	1	-	U	Yes
Hexachloroethane	200	ug/kg	1	-	U	Yes
Indeno(1,2,3-cd)pyrene	40	ug/kg	1	-	U	Yes
Isophorone	80	ug/kg	1	-	U	Yes
1-Methylnaphthalene	80	ug/kg	1	-	U	Yes
2-Methylnaphthalene	80	ug/kg	1	-	U	Yes
2-Nitroaniline	200	ug/kg	1	-	U	Yes
3-Nitroaniline	200	ug/kg	1	-	U	Yes
4-Nitroaniline	200	ug/kg	1	073	U	Yes
Nitrobenzene	80	ug/kg	1	-	U	Yes
N-Nitroso-di-n-propylamine	80	ug/kg	1	-	U	Yes
Nitrosodiphenylamine	200	ug/kg	1	-	υ	Yes
Phenanthrene	40	ug/kg	1	-	U	Yes
Pyrene	40	ug/kg	1	-	U	Yes
1,2,4,5-Tetrachlorobenzene	200	ug/kg	1	2	U	Yes
METHOD: 8	8270D (S	IM)				
Naphthalene	3.9	ug/kg	1	-	U	Yes
1,4-Dioxane	3.9	ug/kg	1	-	U	Yes

Sample location: BMSMC Building 5 Area

Sampling date: 5/26/2016 Matrix: Groundwater

2-Chlorophenol 5.6 ug/l 1 - U Yes 4-Chloro-3-methyl phenol 5.6 ug/l 1 - U Yes 2,4-Dichlorophenol 2.2 ug/l 1 - U Yes 2,4-Dimethylphenol 5.6 ug/l 1 - U Yes 2,4-Dimethylphenol 5.6 ug/l 1 - U Yes 4,6-Dinitro-0-cresol 5.6 ug/l 1 - U Yes 4,6-Dinitro-0-cresol 5.6 ug/l 1 - U Yes 3&4-Methylphenol 2.2 ug/l 1 - U Yes 3&4-Methylphenol 2.2 ug/l 1 - U Yes 3&4-Methylphenol 5.6 ug/l 1 - U Yes 2-Nitrophenol 5.6 ug/l 1 - U Yes 4-Nitrophenol 5.6 ug/l 1 - U Yes 2-Nitrophenol 5.6 ug/l 1 - U Yes 2-Nitrophenol 5.6 ug/l 1 - U Yes 9-Pentachlorophenol 5.6 ug/l 1 - U Yes 9-Pentachloro	Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
2,4-Dichlorophenol 2.2 ug/l 1 - U Yes 2,4-Dimethylphenol 5.6 ug/l 1 - U Yes 2,4-Dinitrophenol 1.1 ug/l 1 - U Yes 4,6-Dinitro-o-cresol 5.6 ug/l 1 - U Yes 2-Methylphenol 2.2 ug/l 1 - U Yes 2-Mitrophenol 5.6 ug/l 1 - U Yes 2-Nitrophenol 5.6 ug/l 1 - U Yes 4-Nitrophenol 5.6 ug/l 1 - U Yes Pentachlorophenol 5.6 ug/l 1 - U Yes Phenol 2.2 ug/l 1 - U Yes 2,4,6-Tetrachlorophenol 5.6 ug/l 1 - U Yes Acenaphthene 1.1 ug/l 1 - U<	•	5.6	ug/l	1	-	U	Yes
2,4-Dimethylphenol 5.6 ug/l 1 - U Yes 2,4-Dinitrophenol 11 ug/l 1 - U Yes 4,6-Dinitro-o-cresol 5.6 ug/l 1 - U Yes 2-Methylphenol 2.2 ug/l 1 - U Yes 2-Nitrophenol 5.6 ug/l 1 - U Yes 4-Nitrophenol 11 ug/l 1 - U Yes 4-Nitrophenol 5.6 ug/l 1 - U Yes 4-Nitrophenol 5.6 ug/l 1 - U Yes Pentachlorophenol 5.6 ug/l 1 - U Yes 2,3,4,6-Tetrachlorophenol 5.6 ug/l 1 - U Yes Acenaphthene 1.1 ug/l 1 - U Yes Acenaphthylene 1.1 ug/l 1 - <td< td=""><td>4-Chloro-3-methyl phenol</td><td>5.6</td><td>ug/l</td><td>1</td><td>-</td><td>U</td><td>Yes</td></td<>	4-Chloro-3-methyl phenol	5.6	ug/l	1	-	U	Yes
2,4-Dinitrophenol 11 ug/l 1 - U Yes 4,6-Dinitro-o-cresol 5.6 ug/l 1 - U Yes 2-Methylphenol 2.2 ug/l 1 - U Yes 3&4-Methylphenol 2.2 ug/l 1 - U Yes 2-Nitrophenol 5.6 ug/l 1 - U Yes 4-Nitrophenol 11 ug/l 1 - U Yes 4-Nitrophenol 5.6 ug/l 1 - U Yes Pentachlorophenol 5.6 ug/l 1 - U Yes Phenol 2.2 ug/l 1 - U Yes 2,4,6-Tetrachlorophenol 5.6 ug/l 1 - U Yes Acenaphthene 1.1 ug/l 1 - U Yes Acenaphthylene 1.1 ug/l 1 - U Yes Acetophenone 2.2 ug/l 1 - U	2,4-Dichlorophenol	2.2	ug/l	1	-	U	Yes
4,6-Dinitro-o-cresol 5.6 ug/l 1 - U Yes 2-Methylphenol 2.2 ug/l 1 - U Yes 3&A-Methylphenol 2.2 ug/l 1 - U Yes 2-Nitrophenol 5.6 ug/l 1 - U Yes 4-Nitrophenol 5.6 ug/l 1 - U Yes Pentachlorophenol 5.6 ug/l 1 - U Yes Phenol 2.2 ug/l 1 - U Yes Phenol 2.2 ug/l 1 - U Yes 2,4,6-Tetrachlorophenol 5.6 ug/l 1 - U Yes Acenaphthene 1.1 ug/l 1 - U Yes Acenaphthylene 1.1 ug/l 1 - U Yes Acetophenone 1.1 ug/l 1 - U Y	2,4-Dimethylphenol	5.6	ug/l	1	-	U	Yes
2-Methylphenol 2.2 ug/l 1 - U Yes 3&4-Methylphenol 2.2 ug/l 1 - U Yes 2-Nitrophenol 5.6 ug/l 1 - U Yes 4-Nitrophenol 5.6 ug/l 1 - U Yes 4-Nitrophenol 5.6 ug/l 1 - U Yes Pentachlorophenol 5.6 ug/l 1 - U Yes Phenol 2.2 ug/l 1 - U Yes 2,3,4,6-Tetrachlorophenol 5.6 ug/l 1 - U Yes 2,3,4,6-Tetrachlorophenol 5.6 ug/l 1 - U Yes 2,4,5-Trichlorophenol 5.6 ug/l 1 - U Yes 2,4,5-Trichlorophenol 5.6 ug/l 1 - U Yes Acenaphthylene 1.1 ug/l 1 - U Yes Acenaphthylene 1.1 ug/l 1 - U Yes Acenaphthylene 1.1 ug/l 1 - U Yes Acetophenone 2.2 ug/l 1 - U Yes Acenaphthylene 1.1 ug/l 1 - U Yes Atrazine 2.2 ug/l 1 - U Yes Benzo(a)anthracene 1.1 ug/l 1 - U Yes Benzo(a)pyrene 1.1 ug/l 1 - U Yes Benzo(a)pyrene 1.1 ug/l 1 - U Yes Benzo(b)fluoranthene 1.1 ug/l 1 - U Yes Benzo(b)fluoranthene 1.1 ug/l 1 - U Yes Benzo(b)fluoranthene 1.1 ug/l 1 - U Yes Benzo(k)fluoranthene	2,4-Dinitrophenol	11	ug/l	1	-	U	Yes
3&4-Methylphenol 2.2 ug/l 1 - U Yes 2-Nitrophenol 5.6 ug/l 1 - U Yes 4-Nitrophenol 11 ug/l 1 - U Yes Pentachlorophenol 5.6 ug/l 1 - U Yes Phenol 2.2 ug/l 1 - U Yes Phenol 5.6 ug/l 1 - U Yes 2.3,4,6-Tetrachlorophenol 5.6 ug/l 1 - U Yes 2.4,5-Trichlorophenol 5.6 ug/l 1 - U Yes Acenaphthene 1.1 ug/l 1 - U Yes Acenaphthene 1.1 ug/l 1 - U Yes Acetophenone 2.2 ug/l 1 - U Yes Acterophenone 1.1 ug/l 1 - U Yes Atrazine 2.2 ug/l 1 - U Yes	4,6-Dinitro-o-cresol	5.6	ug/l	1	-	U	Yes
2-Nitrophenol 5.6 ug/l 1 - U Yes 4-Nitrophenol 11 ug/l 1 - U Yes Pentachlorophenol 5.6 ug/l 1 - U Yes Phenol 2.2 ug/l 1 - U Yes 2,3,4,6-Tetrachlorophenol 5.6 ug/l 1 - U Yes 2,4,5-Trichlorophenol 5.6 ug/l 1 - U Yes 2,4,5-Trichlorophenol 5.6 ug/l 1 - U Yes 2,4,5-Trichlorophenol 5.6 ug/l 1 - U Yes 2,4,6-Trichlorophenol 5.6 ug/l 1 - U Yes Acenaphthene 1.1 ug/l 1 - U Yes Acenaphthene 1.1 ug/l 1 - U Yes Acenaphthylene 1.1 ug/l 1 - U Yes Acetophenone 2.2 ug/l 1 - U Yes Acetophenone 2.2 ug/l 1 - U Yes Atrazine 2.2 ug/l 1 - U Yes Benzaldehyde 5.6 ug/l 1 - U Yes Benzo(a)anthracene 1.1 ug/l 1 - U Yes Benzo(a)pyrene 1.1 ug/l 1 - U Yes Benzo(b)fluoranthene 1.1 ug/l 1 - U Yes Benzo(k)fluoranthene 1.1 ug/l 1 - U Yes Ca-Bonzole 1.1 ug/l 1 - U Yes Carbazole 1.1 ug/	2-Methylphenol	2.2	ug/l	1	-	U	Yes
A-Nitrophenol	3&4-Methylphenol	2.2	ug/l	1	-	U	Yes
Pentachlorophenol 5.6 ug/l 1 - U Yes Phenol 2.2 ug/l 1 - UJ Yes 2,3,4,6-Tetrachlorophenol 5.6 ug/l 1 - U Yes 2,4,5-Trichlorophenol 5.6 ug/l 1 - U Yes 2,4,6-Trichlorophenol 5.6 ug/l 1 - U Yes Acenaphthene 1.1 ug/l 1 - U Yes Acenaphthylene 1.1 ug/l 1 - U Yes Acetophenone 2.2 ug/l 1 - U Yes Acetophenone 2.2 ug/l 1 - U Yes Acetophenone 1.1 ug/l 1 - U Yes Actophonone 1.1 ug/l 1 - U Yes Atrazine 2.2 ug/l 1 - U	2-Nitrophenol	5.6	ug/l	1	-	U	Yes
Phenol 2.2 ug/l 1 - UJ Yes 2,3,4,6-Tetrachlorophenol 5.6 ug/l 1 - U Yes 2,4,5-Trichlorophenol 5.6 ug/l 1 - U Yes 2,4,6-Trichlorophenol 5.6 ug/l 1 - U Yes Acenaphthene 1.1 ug/l 1 - U Yes Acenaphthylene 1.1 ug/l 1 - U Yes Acetophenone 2.2 ug/l 1 - U Yes Acetophenone 2.2 ug/l 1 - U Yes Acetophenone 2.2 ug/l 1 - U Yes Actophenone 2.2 ug/l 1 - U Yes Anthracene 1.1 ug/l 1 - U Yes Benzo(a)anthracene 1.1 ug/l 1 - U	4-Nitrophenol	11	ug/l	1	-	U	Yes
2,3,4,6-Tetrachlorophenol 5.6 ug/l 1 - U Yes 2,4,5-Trichlorophenol 5.6 ug/l 1 - U Yes 2,4,6-Trichlorophenol 5.6 ug/l 1 - U Yes Acenaphthene 1.1 ug/l 1 - U Yes Acenaphthylene 1.1 ug/l 1 - U Yes Acetophenone 2.2 ug/l 1 - U Yes Acetophenone 2.2 ug/l 1 - U Yes Anthracene 1.1 ug/l 1 - U Yes Atrazine 2.2 ug/l 1 - U Yes Benzo(a)purene 1.1 ug/l 1 - U Yes Benzo(a)pyrene 1.1 ug/l 1 - U Yes Benzo(b)fluoranthene 1.1 ug/l 1 - U Yes Benzo(k)fluoranthene 1.1 ug/l 1 -	Pentachlorophenol	5.6	ug/l	1	-	U	Yes
2,4,5-Trichlorophenol 5.6 ug/l 1 - U Yes 2,4,6-Trichlorophenol 5.6 ug/l 1 - U Yes Acenaphthene 1.1 ug/l 1 - U Yes Acenaphthylene 1.1 ug/l 1 - U Yes Acetophenone 2.2 ug/l 1 - U Yes Anthracene 1.1 ug/l 1 - U Yes Antrazine 2.2 ug/l 1 - U Yes Benzaldehyde 5.6 ug/l 1 - U Yes Benzo(a)anthracene 1.1 ug/l 1 - U Yes Benzo(a)pyrene 1.1 ug/l 1 - U Yes Benzo(b)fluoranthene 1.1 ug/l 1 - U Yes Benzo(g,h,i)perylene 1.1 ug/l 1 - U Yes Benzo(k)fluoranthene 1.1 ug/l 1 -	Phenol	2.2	ug/l	1	-	UJ	Yes
2,4,6-Trichlorophenol 5.6 ug/l 1 - U Yes Acenaphthene 1.1 ug/l 1 - U Yes Acenaphthylene 1.1 ug/l 1 - U Yes Acetophenone 2.2 ug/l 1 - U Yes Anthracene 1.1 ug/l 1 - U Yes Atrazine 2.2 ug/l 1 - U Yes Benzaldehyde 5.6 ug/l 1 - U Yes Benzolajanthracene 1.1 ug/l 1 - U Yes Benzolajpyrene 1.1 ug/l 1 - U Yes Benzolajpyrene 1.1 ug/l 1 - U	2,3,4,6-Tetrachlorophenol	5.6	ug/l	1	-	U	Yes
Acenaphthene 1.1 ug/l 1 - U Yes Acenaphthylene 1.1 ug/l 1 - U Yes Acetophenone 2.2 ug/l 1 - U Yes Anthracene 1.1 ug/l 1 - U Yes Atrazine 2.2 ug/l 1 - U Yes Benzaldehyde 5.6 ug/l 1 - U Yes Benzolajanthracene 1.1 ug/l 1 - U <td>2,4,5-Trichlorophenol</td> <td>5.6</td> <td>ug/l</td> <td>1</td> <td>-</td> <td>U</td> <td>Yes</td>	2,4,5-Trichlorophenol	5.6	ug/l	1	-	U	Yes
Acenaphthylene 1.1 ug/l 1 - U Yes Acetophenone 2.2 ug/l 1 - U Yes Anthracene 1.1 ug/l 1 - U Yes Atrazine 2.2 ug/l 1 - U Yes Benzaldehyde 5.6 ug/l 1 - U Yes Benzo(a)anthracene 1.1 ug/l 1 - U Yes Benzo(a)pyrene 1.1 ug/l 1 - U Yes Benzo(b)fluoranthene 1.1 ug/l 1 - U Yes Benzo(k)fluoranthene 1.1 ug/l 1 - U Yes Benzo(k)fluoranthene 1.1 ug/l 1 - U Yes 4-Bromophenyl phenyl ether 1.1 ug/l 1 - U Yes Butyl benzyl phthalate 2.2 ug/l 1 - U Yes 2-Chloroanphthalene 2.2 ug/l 1	2,4,6-Trichlorophenol	5.6	ug/l	1	-	U	Yes
Acenaphthylene 1.1 ug/l 1 - U Yes Acetophenone 2.2 ug/l 1 - U Yes Anthracene 1.1 ug/l 1 - U Yes Atrazine 2.2 ug/l 1 - U Yes Benzaldehyde 5.6 ug/l 1 - U Yes Benzo(a)anthracene 1.1 ug/l 1 - U Yes Benzo(a)pyrene 1.1 ug/l 1 - U Yes Benzo(b)fluoranthene 1.1 ug/l 1 - U Yes Benzo(k)fluoranthene 1.1 ug/l 1 - U Yes Benzo(k)fluoranthene 1.1 ug/l 1 - U Yes 4-Bromophenyl phenyl ether 1.1 ug/l 1 - U Yes Butyl benzyl phthalate 2.2 ug/l 1 - U Yes 2-Chloroanphthalene 2.2 ug/l 1	Acenaphthene	1.1	ug/l	1	-	U	Yes
Anthracene 1.1 ug/l 1 - U Yes Atrazine 2.2 ug/l 1 - U Yes Benzaldehyde 5.6 ug/l 1 - U Yes Benzo(a)anthracene 1.1 ug/l 1 - U Yes Benzo(a)pyrene 1.1 ug/l 1 - U Yes Benzo(b)fluoranthene 1.1 ug/l 1 - U Yes Benzo(k)fluoranthene 1.1 ug/l 1 - U Yes Benzo(k)fluoranthene 1.1 ug/l 1 - U Yes 4-Bromophenyl phenyl ether 1.1 ug/l 1 - U Yes Butyl benzyl phthalate 2.2 ug/l 1 - U Yes 1,1'-Biphenyl 1.1 ug/l 1 - U Yes 2-Chloroaphthalene 2.2 ug/l 1 - U Yes Carbazole 1.1 ug/l 1 - <td>Acenaphthylene</td> <td>1.1</td> <td></td> <td>1</td> <td>-</td> <td>U</td> <td>Yes</td>	Acenaphthylene	1.1		1	-	U	Yes
Anthracene 1.1 ug/l 1 - U Yes Atrazine 2.2 ug/l 1 - U Yes Benzaldehyde 5.6 ug/l 1 - U Yes Benzo(a)anthracene 1.1 ug/l 1 - U Yes Benzo(a)pyrene 1.1 ug/l 1 - U Yes Benzo(b)fluoranthene 1.1 ug/l 1 - U Yes Benzo(k)fluoranthene 1.1 ug/l 1 - U Yes Benzo(k)fluoranthene 1.1 ug/l 1 - U Yes 4-Bromophenyl phenyl ether 1.1 ug/l 1 - U Yes Butyl benzyl phthalate 2.2 ug/l 1 - U Yes 1,1'-Biphenyl 1.1 ug/l 1 - U Yes 2-Chloroaphthalene 2.2 ug/l 1 - U Yes Carbazole 1.1 ug/l 1 - <td>Acetophenone</td> <td>2.2</td> <td>ug/l</td> <td>1</td> <td>-</td> <td>U</td> <td>Yes</td>	Acetophenone	2.2	ug/l	1	-	U	Yes
Benzaldehyde 5.6 ug/l 1 - U Yes Benzo(a)anthracene 1.1 ug/l 1 - U Yes Benzo(a)pyrene 1.1 ug/l 1 - U Yes Benzo(b)fluoranthene 1.1 ug/l 1 - U Yes Benzo(k)fluoranthene 1.1 ug/l 1 - U Yes Benzo(k)fluoranthene 1.1 ug/l 1 - U Yes 4-Bromophenyl phenyl ether 1.1 ug/l 1 - U Yes Butyl benzyl phthalate 2.2 ug/l 1 - U Yes Butyl benzyl phthalate 2.2 ug/l 1 - U Yes 1,1'-Biphenyl 1.1 ug/l 1 - U Yes 2-Chloroaphthalene 2.2 ug/l 1 - U Yes Carbazole 1.1 ug/l 1	Anthracene	1.1		1	-	U	Yes
Benzo(a)anthracene 1.1 ug/l 1 - U Yes Benzo(a)pyrene 1.1 ug/l 1 - U Yes Benzo(b)fluoranthene 1.1 ug/l 1 - U Yes Benzo(k)fluoranthene 1.1 ug/l 1 - U Yes Benzo(k)fluoranthene 1.1 ug/l 1 - U Yes 4-Bromophenyl phenyl ether 1.1 ug/l 1 - U Yes Butyl benzyl phthalate 2.2 ug/l 1 - U Yes Butyl benzyl phthalate 2.2 ug/l 1 - U Yes 1,1'-Biphenyl 1.1 ug/l 1 - U Yes 2-Chloroaphthalene 2.2 ug/l 1 - U Yes Carbazole 1.1 ug/l 1 - U Yes Caprolactam 2.2 ug/l 1	Atrazine	2.2	ug/l	1	-	U	Yes
Benzo(a)pyrene 1.1 ug/l 1 - U Yes Benzo(b)fluoranthene 1.1 ug/l 1 - U Yes Benzo(k)fluoranthene 1.1 ug/l 1 - U Yes Benzo(k)fluoranthene 1.1 ug/l 1 - U Yes 4-Bromophenyl phenyl ether 1.1 ug/l 1 - U Yes Butyl benzyl phthalate 2.2 ug/l 1 - U Yes 1,1'-Biphenyl 1.1 ug/l 1 - U Yes 2-Chloronaphthalene 2.2 ug/l 1 - U Yes 4-Chloroaniline 5.6 ug/l 1 - U Yes Carbazole 1.1 ug/l 1 - U Yes Chrysene 1.1 ug/l 1 - U Yes bis(2-Chloroethoxy)methane 2.2 ug/l 1 - U Yes bis(2-Chloroethyl)ether 2.2 ug/l <	Benzaldehyde	5.6	ug/l	1	-	U	Yes
Benzo(b)fluoranthene 1.1 ug/l 1 - U Yes Benzo(g,h,i)perylene 1.1 ug/l 1 - U Yes Benzo(k)fluoranthene 1.1 ug/l 1 - U Yes 4-Bromophenyl phenyl ether 1.1 ug/l 1 - U Yes Butyl benzyl phthalate 2.2 ug/l 1 - U Yes 1,1'-Biphenyl 1.1 ug/l 1 - U Yes 2-Chloronaphthalene 2.2 ug/l 1 - U Yes 4-Chloroaniline 5.6 ug/l 1 - U Yes Carbazole 1.1 ug/l 1 - U Yes Chrysene 1.1 ug/l 1 - U Yes bis(2-Chloroethoxy)methane 2.2 ug/l 1 - U Yes bis(2-Chloroethyl)ether 2.2 ug/l 1 - U Yes	Benzo(a)anthracene	1.1	ug/l	1	-	U	Yes
Benzo(g,h,i)perylene 1.1 ug/l 1 - U Yes Benzo(k)fluoranthene 1.1 ug/l 1 - U Yes 4-Bromophenyl phenyl ether 1.1 ug/l 1 - U Yes Butyl benzyl phthalate 2.2 ug/l 1 - U Yes 1,1'-Biphenyl 1.1 ug/l 1 - U Yes 2-Chloronaphthalene 2.2 ug/l 1 - U Yes 4-Chloroaniline 5.6 ug/l 1 - U Yes Carbazole 1.1 ug/l 1 - U Yes Caprolactam 2.2 ug/l 1 - U Yes Chrysene 1.1 ug/l 1 - U Yes bis(2-Chloroethoxy)methane 2.2 ug/l 1 - U Yes bis(2-Chloroethyl)ether 2.2 ug/l 1 - UJ Yes	Benzo(a)pyrene	1.1	ug/l	1	-	U	Yes
Benzo(k)fluoranthene 1.1 ug/l 1 - U Yes 4-Bromophenyl phenyl ether 1.1 ug/l 1 - U Yes Butyl benzyl phthalate 2.2 ug/l 1 - U Yes 1,1'-Biphenyl 1.1 ug/l 1 - U Yes 2-Chloronaphthalene 2.2 ug/l 1 - U Yes 4-Chloroaniline 5.6 ug/l 1 - U Yes Carbazole 1.1 ug/l 1 - U Yes Caprolactam 2.2 ug/l 1 - U Yes Chrysene 1.1 ug/l 1 - U Yes bis(2-Chloroethoxy)methane 2.2 ug/l 1 - U Yes bis(2-Chloroethyl)ether 2.2 ug/l 1 - U Yes	Benzo(b)fluoranthene	1.1	ug/l	1	-	U	Yes
4-Bromophenyl phenyl ether 1.1 ug/l 1 - U Yes Butyl benzyl phthalate 2.2 ug/l 1 - U Yes 1,1'-Biphenyl 1.1 ug/l 1 - U Yes 2-Chloronaphthalene 2.2 ug/l 1 - U Yes 4-Chloroaniline 5.6 ug/l 1 - U Yes Carbazole 1.1 ug/l 1 - U Yes Caprolactam 2.2 ug/l 1 - U Yes Chrysene 1.1 ug/l 1 - U Yes bis(2-Chloroethoxy)methane 2.2 ug/l 1 - U Yes bis(2-Chloroethyl)ether 2.2 ug/l 1 - UJ Yes	Benzo(g,h,i)perylene	1.1	ug/l	1	-	U	Yes
Butyl benzyl phthalate 2.2 ug/l 1 - U Yes 1,1'-Biphenyl 1.1 ug/l 1 - U Yes 2-Chloronaphthalene 2.2 ug/l 1 - U Yes 4-Chloroaniline 5.6 ug/l 1 - U Yes Carbazole 1.1 ug/l 1 - U Yes Caprolactam 2.2 ug/l 1 - U Yes Chrysene 1.1 ug/l 1 - U Yes bis(2-Chloroethoxy)methane 2.2 ug/l 1 - U Yes bis(2-Chloroethyl)ether 2.2 ug/l 1 - UJ Yes	Benzo(k)fluoranthene	1.1	ug/l	1	-	U	Yes
1,1'-Biphenyl 1.1 ug/l 1 - U Yes 2-Chloronaphthalene 2.2 ug/l 1 - U Yes 4-Chloroaniline 5.6 ug/l 1 - U Yes Carbazole 1.1 ug/l 1 - U Yes Caprolactam 2.2 ug/l 1 - U Yes Chrysene 1.1 ug/l 1 - U Yes bis(2-Chloroethoxy)methane 2.2 ug/l 1 - U Yes bis(2-Chloroethyl)ether 2.2 ug/l 1 - UJ Yes	4-Bromophenyl phenyl ether	1.1	ug/l	1	-	U	Yes
2-Chloronaphthalene 2.2 ug/l 1 - U Yes 4-Chloroaniline 5.6 ug/l 1 - U Yes Carbazole 1.1 ug/l 1 - U Yes Caprolactam 2.2 ug/l 1 - U Yes Chrysene 1.1 ug/l 1 - U Yes bis(2-Chloroethoxy)methane 2.2 ug/l 1 - UJ Yes bis(2-Chloroethyl)ether 2.2 ug/l 1 - UJ Yes	Butyl benzyl phthalate	2.2	ug/l	1	-	U	Yes
4-Chloroaniline 5.6 ug/l 1 - U Yes Carbazole 1.1 ug/l 1 - U Yes Caprolactam 2.2 ug/l 1 - U Yes Chrysene 1.1 ug/l 1 - U Yes bis(2-Chloroethoxy)methane 2.2 ug/l 1 - U Yes bis(2-Chloroethyl)ether 2.2 ug/l 1 - UJ Yes	1,1'-Biphenyl	1.1	ug/l	1	-	U	Yes
Carbazole 1.1 ug/l 1 - U Yes Caprolactam 2.2 ug/l 1 - U Yes Chrysene 1.1 ug/l 1 - U Yes bis(2-Chloroethoxy)methane 2.2 ug/l 1 - U Yes bis(2-Chloroethyl)ether 2.2 ug/l 1 - UJ Yes	2-Chloronaphthalene	2.2	ug/l	1	-	U	Yes
Caprolactam 2.2 ug/l 1 - U Yes Chrysene 1.1 ug/l 1 - U Yes bis(2-Chloroethoxy)methane 2.2 ug/l 1 - U Yes bis(2-Chloroethyl)ether 2.2 ug/l 1 - UJ Yes	4-Chloroaniline	5.6	ug/l	1	-	U	Yes
Chrysene 1.1 ug/l 1 - U Yes bis(2-Chloroethoxy)methane 2.2 ug/l 1 - U Yes bis(2-Chloroethyl)ether 2.2 ug/l 1 - UJ Yes	Carbazole	1.1	ug/l	1	-	U	Yes
bis(2-Chloroethoxy)methane 2.2 ug/l 1 - U Yes bis(2-Chloroethyl)ether 2.2 ug/l 1 - UJ Yes	Caprolactam	2.2	ug/l	1	-	U	Yes
bis(2-Chloroethyl)ether 2.2 ug/l 1 - UJ Yes	Chrysene	1.1	ug/l	1	-	Ų	Yes
	bis(2-Chloroethoxy)methane	2.2	ug/l	1	-	U	Yes
bis(2-Chloroisopropyl)ether 2.2 ug/l 1 - U Yes	bis(2-Chloroethyl)ether	2.2	ug/l	1	-	UJ	Yes
	bis(2-Chloroisopropyl)ether	2.2	ug/l	1	-	U	Yes

4-Chlorophenyl phenyl ether	2.2	ug/l	1	-	U	Yes		
2,4-Dinitrotoluene	1.1	ug/l	1	-	U	Yes		
2,6-Dinitrotoluene	1.1	ug/l	1	-	U	Yes		
3,3'-Dichlorobenzidine	2.2	ug/l	1		U	Yes		
Dibenzo(a,h)anthracene	1.1	ug/l	1	2	U	Yes		
Dibenzofuran	5.6	ug/l	1		U	Yes		
Di-n-butyl phthalate	2.2	ug/l	1	-	U	Yes		
Di-n-octyl phthalate	2.2	ug/l	1		U	Yes		
Diethyl phthalate	2.2	ug/l	1		U	Yes		
Dimethyl phthalate	2.2	ug/l	1		Ų	Yes		
bis(2-Ethylhexyl)phthalate	2.2	ug/l	1	-	U	Yes		
Fluoranthene	1.1	ug/l	1		U	Yes		
Fluorene	1.1	ug/l	1	-	Ų	Yes		
Hexachlorobenzene	1.1	ug/l	1	-	U	Yes		
Hexachlorobutadiene	1.1	ug/l	1	-	U	Yes		
Hexachlorocyclopentadiene	11	ug/l	1		U	Yes		
Hexachloroethane	2.2	ug/l	1	-	U	Yes		
Indeno(1,2,3-cd)pyrene	1.1	ug/l	1	-	U	Yes		
Isophorone	2.2	ug/l	1	2	U	Yes		
1-Methylnaphthalene	1.1	ug/l	1	-	U	Yes		
2-Methylnaphthalene	1.1	ug/l	1	17.1	U	Yes		
2-Nitroaniline	5.6	ug/l	1	-	UJ	Yes		
3-Nitroaniline	5.6	ug/l	1	-	UJ	Yes		
4-Nitroaniline	5.6	ug/l	1	-	U	Yes		
Nitrobenzene	2.2	ug/l	1	-	U	Yes		
N-Nitroso-di-n-propylamine	2.2	ug/l	1	-	UJ	Yes		
Nitrosodiphenylamine	5.6	ug/l	1	-	U	Yes		
Phenanthrene	1.1	ug/l	1	-	U	Yes		
Pyrene	1.1	ug/l	1	25/4	U	Yes		
1,2,4,5-Tetrachlorobenzene	2.2	ug/l	1		U	Yes		
METHOD: 8270D (SIM)								
Naphthalene	0.11	ug/l	1	121	U	Yes		
1,4-Dioxane	0.298	ug/i	1		Ų.	Yes		
wy i wrociatio	0.230	ug/1	-	-	- 5	1.62		

.

Sample location: BMSMC Building 5 Area

Sampling date: 5/26/2016 Matrix: Groundwater

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
2-Chlorophenol	5.6	ug/l	1	-	U	Yes
4-Chloro-3-methyl phenol	5.6	ug/l	1	-	U	Yes
2,4-Dichlorophenol	2.2	ug/l	1	-	U	Yes
2,4-Dimethylphenol	5.6	ug/l	1	-	U	Yes
2,4-Dinitrophenol	11	ug/l	1	-	U	Yes
4,6-Dinitro-o-cresol	5.6	ug/l	1	-	U	Yes
2-Methylphenol	2.2	ug/l	1	-	U	Yes
3&4-Methylphenol	2.2	ug/l	1	-	U	Yes
2-Nitrophenol	5.6	ug/l	1	-	U	Yes
4-Nitrophenol	11	ug/l	1	-	U	Yes
Pentachlorophenol	5.6	ug/l	1	-	U	Yes
Phenol	2.2	ug/l	1	4	IJ	Yes
2,3,4,6-Tetrachlorophenol	5.6	ug/l	1	-	U	Yes
2,4,5-Trichlorophenol	5.6	ug/l	1	-	U	Yes
2,4,6-Trichlorophenol	5.6	ug/l	1	-	U	Yes
Acenaphthene	1.1	ug/l	1	-	U	Yes
Acenaphthylene	1.1	ug/l	1	-	U	Yes
Acetophenone	2.2	ug/l	1	-	U	Yes
Anthracene	1.1	ug/l	1	-	U	Yes
Atrazine	2.2	ug/l	1	-	U	Yes
Benzaldehyde	5.6	ug/l	1	-	U	Yes
Benzo(a)anthracene	1.1	ug/l	1	-	U	Yes
Benzo(a)pyrene	1.1	ug/l	1	-	U	Yes
Benzo(b)fluoranthene	1.1	ug/l	1	-	U	Yes
Benzo(g,h,i)perylene	1.1	ug/l	1	-	U	Yes
Benzo(k)fluoranthene	1.1	ug/l	1	-	U	Yes
4-Bromophenyl phenyl ether	1.1	ug/l	1	-	U	Yes
Butyl benzyl phthalate	2.2	ug/l	1	-	U	Yes
1,1'-Biphenyl	1.1	ug/l	1	-	U	Yes
2-Chloronaphthalene	2.2	ug/l	1	-	U	Yes
4-Chloroaniline	5.6	ug/l	1	-	U	Yes
Carbazole	1.1	ug/l	1	-	U	Yes
Caprolactam	2.2	ug/l	1	-	U	Yes
Chrysene	1.1	ug/l	1	-	U	Yes
bis(2-Chloroethoxy)methane	2.2	ug/l	1	-	U	Yes
bis(2-Chloroethyl)ether	2.2	ug/l	1	-	UJ	Yes
bis(2-Chloroisopropyl)ether	2.2	ug/l	1	-	U	Yes
4-Chlorophenyl phenyl ether	2.2	ug/l	1	-	U	Yes

2,4-Dinitrotoluene	1.1	ug/l	1	-	U	Yes
2,6-Dinitrotoluene	1.1	ug/l	1	-	U	Yes
3,3'-Dichlorobenzidine	2.2	ug/l	1	-	U	Yes
Dibenzo(a,h)anthracene	1.1	ug/l	1	-	U	Yes
Dibenzofuran	5.6	ug/l	1	-	U	Yes
Di-n-butyl phthalate	2.2	ug/l	1	-	U	Yes
Di-n-octyl phthalate	2.2	ug/l	1		U	Yes
Diethyl phthalate	2.2	ug/l	1	-	U	Yes
Dimethyl phthalate	2.2	ug/i	1	-	U	Yes
bis(2-Ethylhexyl)phthalate	2.2	ug/l	1	(5)	U	Yes
Fluoranthene	1.8	ug/l	1	-	-	Yes
Fluorene	1.1	ug/l	1	-	U	Yes
Hexachlorobenzene	1.1	ug/l	1	-	U	Yes
Hexachlorobutadiene	1.1	ug/l	1	-	U	Yes
Hexachlorocyclopentadiene	11	ug/l	1	-	U	Yes
Hexachloroethane	2.2	ug/l	1	2	U	Yes
Indeno(1,2,3-cd)pyrene	1.1	ug/l	1	-	U	Yes
Isophorone	2.2	ug/l	1		U	Yes
1-Methylnaphthalene	1.1	ug/l	1	-	U	Yes
2-Methylnaphthalene	1.1	ug/l	1	-	U	Yes
2-Nitroaniline	5.6	ug/l	1	-	נט	Yes
3-Nitroaniline	5.6	ug/l	1	-	IJ	Yes
4-Nitroaniline	5.6	ug/l	1	-	U	Yes
Nitrobenzene	2.2	ug/l	1	-	U	Yes
N-Nitroso-di-n-propylamine	2.2	ug/l	1	-	UJ	Yes
Nitrosodiphenylamine	5.6	ug/l	1	376	U	Yes
Phenanthrene	1.1	ug/l	1	-	U	Yes
Pyrene	1.3	ug/l	1	-	-	Yes
1,2,4,5-Tetrachlorobenzene	2.2	ug/l	1	-	U	Yes
METHOD:	8270D (SI	M)				
Naphthalene	0.11	ug/l	1	2	U	Yes
1,4-Dioxane	0.11	ug/l	1	-	U	Yes

* 1 2

[Project Number:_JC21103 Date:May_25-26,_2016 Shipping Date:May_26,_2016 EPA Region:2
REVIEW OF SEMIVOLATILE OF	-
The following guidelines for evaluating volatile organization actions. This document will assist the remake more informed decision and in better serving results were assessed according to USEPA data following order of precedence: EPA Hazardous W. 2015 –Revision 0. Semivolatile Data Validation. The QC on the data review worksheets are from the primar noted.	viewer in using professional judgment to the needs of the data users. The sample validation guidance documents in the aste Support Section, SOP HW-35A, July criteria and data validation actions listed
The hardcopied (laboratory name) _Accutest	data package received has been summarized. The data review for SVOCs
Lab. Project/SDG No.:JC21103 No. of Samples:8_Full_scan/8_SIM	Sample matrix: _Soil/Groundwater
Trip blank No.: - Field blank No.: - JC21103-2 Field duplicate No.: JC21103-3/JC21103-4	
X Data CompletenessX Holding TimesX GC/MS TuningX Internal Standard PerformanceX BlanksX Surrogate RecoveriesX Matrix Spike/Matrix Spike Duplicate	X Laboratory Control SpikesX Field DuplicatesX CalibrationsX Compound IdentificationsX Compound QuantitationX Quantitation Limits
Overall Comments:_ABN_TCL_list_by_method_SW846-82 _analyzed_by_method_SW846-8270D_(SIM)	!70D;_Naphthalene_and_1,4-Dioxane_
Definition of Qualifiers:	
J- Estimated results U- Compound not detected R- Rejected data UJ- Estimated nondetect Reviewer:	

DATA COMPLETENESS

MISSING INFORMATION	DATE LAB. CONTACTED	DATE RECEIVED

All criteria were met _X	
Criteria were not met	
and/or see below	

HOLDING TIMES

The objective of this parameter is to ascertain the validity of the results based on the holding time of the sample from time of collection to the time of analysis.

Complete table for all samples and note the analysis and/or preservation not within criteria

SAMPLE ID	DATE SAMPLED	DATE EXTRACTED/ANALYZED	рН	ACTION
All samples extract	ed and analyzed w	ithin method recommended ho	olding	time. Sample preservation was acceptable.
Cooler temperat	ture (Criteria: 4	+2°C)· 4°C		

<u>Actions</u>

Results will be qualified based on the criteria of the following Table:

Table 1. Holding Time Actions for Semivolatile Analyses

		Ing Time Actions for Semic		tion	
Matrix	Preserved	Criteria	Detected Associated Compounds	Non-Detected Associated Compounds	
No		≤7 days (for extraction) ≤40 days (for analysis)	Use professi	onal judgment	
	No	> 7 days (for extraction) > 40 days (for analysis)	.1	Use professional judgment	
Aqueous	Yes	≤ 7 days (for extraction) ≤ 40 days (for analysis)	No qualification		
	Yes	> 7 days (for extraction) > 40 days (for analysis)	1.	ÜJ	
	Yes/No	Grossly Exceeded	J	UJ or R	
	No	≤ 14 days (for extraction) ≤ 40 days (for analysis)	Use profession	onal judgment	
Non-Aqueous	No	> 14 days (for extraction) > 40 days (for analysis)	J Us profess juden		
	Yes	≤ 14 days (for extraction) ≤ 40 days (for analysis)	No qualification		
	Yes	> 14 days (for extraction) > 40 days (for analysis)	J	ÜJ	
	Yes/No	Grossly Exceeded	Į.	UJ or R	

All criteria were metX Criteria were not met see below
GC/MS TUNING
The assessment of the tuning results is to determine if the sample instrumentation is within the standard tuning QC limits
_X The DFTPP performance results were reviewed and found to be within the specified criteria.
_X DFTPP tuning was performed for every 12 hours of sample analysis.
If no, use professional judgment to determine whether the associated data should be accepted, qualified or rejected.
Notes: These requirements do not apply when samples are analyzed by the Selected Ion Monitoring (SIM) technique.
All mass spectrometer conditions must be identical to those used during the sample analysis. Background subtraction actions resulting in spectral distortion are unacceptable Notes: No data should be qualified based of DFTPP failure.
The requirement to analyze the instrument performance check solution is optional when analysis of PAHs/pentachlorophenol is to be performed by the SIM technique.

LIST	tne	samples	affected:

Actions:

- 1. If sample are analyzed without a preceding valid instrument performance check or are analyzed 12 hours after the Instrument Performance Check, qualify all data in those samples as unusable
- 2. If ion abundance criteria are not met, use professional judgment to determine to what extent the data may be utilized.
- 3. State in the Data Review Narrative, decisions to use analytical data associated with DFTPP instrument performance checks not meeting the contract requirements.
- 4. Use professional judgment to determine if associated data should be qualified based on the spectrum of the mass calibration compounds.

All criteria were met	_X_	
Criteria were not met		
and/or see below		

INITIAL CALIBRATION VERIFICATION

Compliance requirements for satisfactory instrument calibration are established to ensure that the instrument is capable of producing and maintaining acceptable quantitative data.

Date of initial calibration:05/24-25/16_(Scan)_	05/17/16_(SIM)		
Instrument ID numbers:GCMS3E	GCMS3M		
Matrix/Level:Aqueous/low	Aqueous/low		
Date of initial calibration:05/13/2016_(Scan) Instrument ID numbers:GCMS6P	05/31/16;06/06/16(SIM) GCMS4P		
Matrix/Level:Aqueous/low	Aqueous/low		
Date of initial calibration:_05/17-18//2016_(Scan)_ Instrument ID numbers:GCMSZ	06/04/16;06/05/16(Scan) GCMS3P		
Matrix/Level:Aqueous/low			
DATE LAB FILE CRITERIA OUT COMPOL	IND SAMPLES		

	ATE	LAB ID#	FILE	CRITERIA OUT RFs, %RSD, %D, r	COMPOUND	SAMPLES AFFECTED	
-	Initial a	and initi	al calib	ration verification mee	ts the method and guida	ance validation document	
	performance criteria.						

Actions:

Qualify the initial calibration analytes listed in Table 2 using the following criteria:

Table 3. Initial Calibration Actions for Semivolatile Analysis

Criteria	Action		
Criteria	Detect	Non-detect	
Initial Calibration not performed at specified frequency and sequence	Use professional judgment R	Use professional judgment R	
Initial Calibration not performed at the specified concentrations	Į.	ÜJ	
RRF < Minimum RRF in Table 2 for target analyte	Use professional judgment J L or R	R	
RRF > Minimum RRF in Table 2 for target analyte	No qualification	No qualification	
%RSD > Maximum %RSD in Table 2 for target analyte	J	Use professional judgment	
%RSD ≤ Maximum %RSD in Table 2 for target analyte	No qualification	No qualification	

Initial Calibration

Table 2. RRF, %RSD, and %D Acceptance Criteria in Initial Calibration and CCV for Semivolatile Analysis

Analyte	Minimum RRF	Maximum %RSD	Opening Maximum %D ¹	Opening Maximum %D ^t
1,4-Dioxane	0.010	40.0	-40.0	- 50.0
Benzaldehyde	0.100	40.0	-40.0	- 50.0
Phenol	0.080	20,0	-20.0	- 25.0
Bis(2-chloroethyl)ether	0.100	20.0	-20.0	-25.0
2-Chlorophenol	0.200	20.0	-20.0	- 25.0
2-Methylphenol	0.010	20.0	= 20.0	- 25.0
3-Methylphenol	0.010	20.0	-20.0	- 25.0
2,2'-Oxybis-(1-chloropropane)	0.010	20.0	€25.0	-50.0
Acetophenone	0.060	20,0	-20.0	-25.0
4-Methylphenol	0,010	20.0	-20.0	-25.0
N-Nitroso-di-n-propylamine	0.080	20.0	- 25.0	-25.0
llevachloroethane	0.100	20,0	-20.0	-25.0
Nitrobenzene	0.090	20.0	-20.0	-25.0
Sophorone	0.100	20.0	- 20.0	-25.0
2-Nitrophenol	0.060	20.0	-20.0	-25.0
2,4-Dimethylphenol	0.050	20.0	-25.0	-50.0
Bis(2-chloroethoxy)methane	0.080	20.0	-20,0	-25.0
2,4-Dichlorophenol	0.060	20.0	=20.0	-25.0
Naphthalene	0.200	20.0	-20.0	-25.0
1-Chloroaniline	0.010	40.0	-40.0	-50.0
Texachlorobutadiene	0.040	20.0	-20.0	- 25.0
Caprolactam	0.010	40.0	= 30.0	- 50.0
4-Chloro-3-methylphenol	0.040	20.0	-20.0	= 25.0
2-Methylnaphthalene	0.100	20.0	-20.0	-25.0
lexachlorocyclopentadiene	0.010	40.0	-40.0	-50.0
2,4,6-Trichlorophenol	0.090	20.0	-20.0	-25.0
2,4.5-Trichlorophenol	0.100	20.0	-20.0	-25.0
1,1'-Biphenyl	0,200	20,0	-20.0	- 25.0

Analyte	Minimum RRF	Maximum %RSD	Opening Maximum %D ¹	Opening Maximum %D [‡]
2-Chloronaphthalene	0.300	20.0	- 20.0	-25.0
2-Nitroaniline	0.060	20.0	- 25.0	- 25.0
Dimethy Iphthalate	0.300	20.0	- 25.0	-25.0
2.6-Dinitrotoluene	0.080	20.0	-20.0	- 25.0
Acenaphthylene	0.400	20.0	-20.0	-25.0
3-Nitroaniline	0.010	20.0	-25.0	-50.0
Acenaphthene	0,200	20.0	-20.0	-25.0
2.4-Dinitrophenol	0.010	40.0	- 50.0	- 50.0
4-Nitrophenol	0.010	40,0	-40.0	- 50.0
Dibenzofuran	0.300	20.0	~ 20.0	- 25.0
2,4-Dinitrotoluene	0.070	20.0	-20.0	-25.0
Diethylphthalate	0.300	20.0	-20.0	= 25.0
1.2,4,5-Tetrachlorobenzene	0.100	20.0	-20.0	-25.0
4-Chlorophenyl-phenylether	0.100	20,0	- 20.0	-25.0
Fluorene	0.200	20.0	-20.0	- 25.0
4-Nitroaniline	0.010	40,0	-40.0	-50.0
4,6-Dinitro-2-methylphenol	0.010	40,0	-30.0	-50.0
4-Bromophenyl-phenyl ether	0.070	20.0	-20.0	= 25.0
N-Nitrosodiphenylamine	0.100	20.0	±20.0	-25.0
Hexachlorobenzene	0.050	20.0	-20.0	-25.0
Atrazine	0.010	40,0	-25.0	- 50.0
Pentachlorophenol	0.010	40.0	-40.0	- 50.0
Phenanthrene	0.200	20.0	-20.0	-25.0
Anthracene	0.200	20.0	-20.0	-25.0
Carbazole	0.050	20.0	-20.0	-25.0
Di-n-butylphthalate	0.500	20.0	-20.0	-25.0
Fluoranthene	0.100	20,0	- 20,0	- 25.0
Pyrene	0,400	20.0	-25.0	-50.0
Butylbenzylphthalate	0.100	20.0	-25.0	- 50.0

Analyte	Minimum RRF	Maximum %RSD	Opening Maximum %D ¹	Opening Maximum %D ¹
3,3'-Dichlorobenzidine	0.010	40,0	-40.0	-50.0
Benzo(a)anthracene	0.300	20,0	= 20.0	= 25.0
Chrysene	0.200	20.0	- 20.0	- 50.0
Bis(2-ethylhexyl) phthalate	0.200	20.0	-25.0	-50.0
Di-n-octylphthalate	0.010	40,0	-40.0	-50.0
Benzo(b)fluoranthene	0.010	20.0	-25.0	- 50,0
Benzo(k)fluoranthene	0.010	20,0	-25.0	-50.0
Benzo(a)pyrene	0.010	20.0	- 20.0	50.0
Indeno(1,2,3-cd)pyrene	0.010	20.0	= 25.0	-50.0
Dibenzo(a,h)anthracene	0.010	20.0	= 25.0	÷ 50.0
Benzo(g,h,i)perylene	0.010	20.0	- 30.0	- 50.0
2,3,4,6-Tetrachlorophenol	0.040	20.0	= 20.0	= 50.0
Naphthalene	0.600	20.0	= 25.0	- 25.0
2-Methylnaphthalene	0.300	20.0	- 20.0	-25.0
Acenaphthylene	0.900	20.0	= 20.0	- 25.0
Acenaphthene	0.500	20,0	-20.0	- 25.0
Fluorene	0.700	20.0	± 25.0	50.0
Phenanthrene	0.300	20.0	= 25.0	- 50.0
Anthracene	0.400	20.0	= 25.0	= 50.0
Fluoranthene	0.400	20.0	= 25.0	- 50.0
Pyrene	0.500	:20.0	=30.0	-50.0
Benzo(a)anthracene	0,400	20.0	= 25.0	- 50.0
Chyrsene	0.400	20.0	= 25.0	-50.0
Benzo(b)fluoranthene	0,100	20.0	-30.0	- 50.0
Benzo(k)fluoranthene	0.100	20.0	-30.0	= 50.0
Benzo(a)pyrene	0.100	20.0	-25.0	= 50.0
ndeno(1,2,3-ed)pyrene	0,100	20.0	=40.0	- 50.0
Dibenzo(a,h)anthracene	0.010	25.0	- 40.0	- 50.0
Benzo(g,h,i)perylene	0.020	25.0	- 40.0	- 50.0

Pentachlorophenol	0.010	40.0	50,0	-50.0	٦
Deuterated Monitoring Compounds					

Analyte	Minimum RRF	Maximum %RSD	Opening Maximum %D ¹	Closing Maximum %D
1,4-Dioxane-d ₈	0.010	20.0	= 25.0	-50.0
Phenol-ds	0.010	20,0	-25.0	=25.0
Bis-(2-chloroethyl)ether-d _s	0.100	20.0	= 20.0	-25.0
2-Chlorophenol-di	0.200	20.0	-20.0	-25.0
1-Methylphenol-ds	0.010	20.0	-20,0	-25.0
4-Chloroaniline-d ₁	0.010	40.0	-40.0	-50.0
Nitrobenzene-d ₅	0.050	20.0	= 20.0	- 25.0
2-Nitrophenol-d ₄	0.050	20.0	-20.0	-25.0
2,4-Dichlorophenol-d;	0.060	20.0	= 20,0	-25.0
Dimethylphthalate-d ₆	0,300	20,0	-20.0	-25.0
Acenaphthylene-d _x	0,400	20.0	+ 20.0	-25.0
4-Nitrophenol-d ₄	0.010	40.0	-40.0	-50.0
Fluorene-d _{in}	0,100	20.0	= 20.0	± 25.0
4,6-Dinitro-2-methylphenol-dy	0.010	40.0	=30,0	-50.0
Anthracene-d ₁₀	0.300	20.0	- 20.0	-25.0
Pyrene-d ₁₀	0.300	20.0	-25.0	- 50.0
Benzo(a)pyrene-d ₁₂	0.010	20.0	-20.0	- 50.0
Fluoranthene-d _{in} (SIM)	0.400	20.0	±25.0	- 50.0
2-Methylnaphthalene-d ₁₀ (SIM)	0.300	20.0	- 20.0	- 25.0

If a closing CCV is acting as an opening CCV, all target analytes must meet the requirements for an opening CCV.

Note: If analysis by SIM technique is requested for PAH/pentachlorophenols, calibration standards analyzed at 0.10, 0.20, 0.40, 0.80, and 1.0 ng/uL for each target compound of interest and the associated DMCs. Pentachlorophenol will require only a four point initial calibration at 0.20, 0.40, 0.80, and 1.0 ng/uL.

All criteria were met
Criteria were not met
and/or see belowX

CONTINUING CALIBRATION VERIFICATION

Compliance requirements for satisfactory instrument calibration are established to ensure that the instrument is capable of producing and maintaining acceptable quantitative data.

Date of initial calibration:05/24-25/16_(Scan)	05/17/16(SIM)
Date of initial calibration verification (ICV):_05/25/16	05/17-18/16
Date of continuing calibration verification (CCV):_06/03/16	06/01/16
06/10/16	
Date of closing CCV:	-
Instrument ID numbers:GCMS3E	GCMS3M
Matrix/Level:Aqueous/low	Aqueous/low
Date of initial calibration:06/04-05/16_(Scan)	05/31/16;_06/06/16_(SIM)
Date of initial calibration verification (ICV):06/05-06/16	06/01/16;_06/06/16
Date of continuing calibration verification (CCV):06/13/16	06/01/16;_06/03/16;_06/06-07/16_
Date of closing CCV:	-
Instrument ID numbers: GCMS3P	GCMS4P
Matrix/Level:Aqueous/low	
Date of initial calibration:5/17-18/16_(Scan)	05/13/16_(Scan)
Date of initial calibration verification (ICV):05/18/16	05/13/16;_05/16/16
Date of continuing calibration verification (CCV):05/27/16	
Date of closing CCV:	-
Instrument ID numbers:GCMSZ	GCMS6P
Matrix/Level:Aqueous/low	

DATE	LAB FILE	CRITERIA OUT	COMPOUND	SAMPLES
DAIL	ID#	RFs, %RSD, <u>%D</u> , r	COMITOUND	AFFECTED
GCMS6P				3
06/01/16	cc1209-50	22.1	Hexachlorocyclopentadiene*	JC21103-4; -6
		-27.4	2-nitroaniline	
		-22.4	Butylbenzylphthalate	
		-32.3	di-n-octylphthalate*	
GCM3E				
06/10/16	cc3662-50	44.3	4-nitrophenol	JC21103-3
		29.0 Hexachlorocyclopentadiene		
		43.9	Pentachlorophenol	
GCMS3P				
06/13/16	CC2460-50	-32.0	2-nitrophenol	JC21103-5
	:	-20.8	Hexachlorobutadiene	
	į	-29.7	2-nitroaniline	
		-30.6	2,6-dinitrotoluene	
		-22.1	2-nitroaniline	
		-31.8	2,4-dinitrophenol*	

DATE	LAB FILE ID#	CRITERIA OUT RFs, %RSD, % D , r	COMPOUND	SAMPLES AFFECTED
GCMS3P	IDIF	111 3, 701 00, 700, 1		ATTLOTED
06/13/16	CC2460-50	-23.5	4-nitrophenol*	JC21103-5
		-24.8	2,3,4,6-pentachlorophenol	
		-20.3	4,6-dinitro-o-cresol*	
		-21.6	3,3'-dichlorobenzidine*	
GCMS4P				
06/01/16	icv879-1.0	-32.8	1,4-dioxane*	JC21103-1; -2; -8
06/03/16	cc879-1.0	-23.1	1,4-dioxane*	JC21103-3; -4; -5; -
				6
GCMSZ				
05/28/16	cc5533-50	-29.2	phenol	JC21103-1; -2; -7;
		-24.7	Bis(2-chloroethyl)ether	-8
		-25.9	N-nitroso-di-n-propylamine	
		-39.5	Caprolactam*]
		-28.9	2-nitroaniline]
		-22.9	3-nitroaniline]
		-24.8	2,4-dinitrophenol*]
		-21.8	4,6-dinitro-o-cresol*]

Note: Initial and continuing calibration verifications meet the method and guidance document required performance criteria except the cases describe in this document.

No closing calibration verification included in data package. No action taken, professional judgment.

QC samples are not validated.

Actions:

Notes: Verify that the CCV is run at the required frequency (an opening and closing CCV must be run within 12-hour period).

All DMCs must meet the RRF values given in Table 2. No qualification of the data is necessary on DMCs RRF and %RSD/%D alone. Use professional judgment to evaluate DMCs and %RSD/%D data in conjunction with DMCs recoveries to determine the need for qualification of the data.

Qualify the initial calibration analytes listed in Table 2 using the following criteria in the CCVs:

^{*} Analytes with % difference in the continue calibration verification outside the method performance criteria but within the validation guidelines criteria, +40 %. No action taken.

Table 4. CCV Actions for Semivolatile Analysis

Critaria fan Onanina CCV	Cuitoria for Clavina CCV	Action		
Criteria for Opening CCV	Criteria for Closing CCV -	Detect	Non-detect	
CCV not performed at required frequency and sequence	CCV not performed at required frequency	Use professional judgment R	Use professional judgment R	
CCV not performed at specified concentration	CCV not performed at specified concentration	Use professional judgment	Use professional judgment	
RRF < Minimum RRF in Table 2 for target analyte	RRF < Minimum RRF in Table 2 for target analyte	Use professional judgment J or R	R	
RRF > Minimum RRF in Table 2 RRF > Minimum RRF in Table 2 for target analyte for target analyte		No qualification	No qualification	
%D outside the Opening Maximum %D limits in Table 2 for target analyte %D outside the Closing Maximum %D limits in Table 2 for target analyte		J	(j	
O within the inclusive Opening iximum %D limits in Table 2 Maximum %D limits in Table 2 for target analyte		No qualification	No qualification	

All criteria were metX	
Criteria were not met	
and/or see below	

BLANK ANALYSIS RESULTS (Sections 1 & 2)

The assessment of the blank analysis results is to determine the existence and magnitude of contamination problems. The criteria for evaluation of blanks apply only to blanks associated with the samples, including trip, equipment, and laboratory blanks. If problems with any blanks exist, all data associated with the case must be carefully evaluated to determine whether or not there is an inherent variability in the data for the case, or if the problem is an isolated occurrence not affecting other data.

List the contamination in the blanks below. High and low levels blanks must be treated separately.

Notes: The concentration of non-target compounds in all blanks must be less than or equal to 10 ug/L.

The concentration of target compounds in all blanks must be less than its CRQL listed in the method.

Samples taken from a drinking water tap do not have and associated field blank.

Laboratory blanks

DATE ANALYZED	LAB ID	LEVEL/ MATRIX	COMPOUND	CONCENTRATION UNITS
	3.22	80 323		
Field/ <u>Equipmen</u>				
DATE Analyzed	LABID	LEVEL/ MATRIX	COMPOUND	CONCENTRATION UNITS
_data_package.			nent_blankNo_field/trip 	
			**	

All criteria were met _X_	
Criteria were not met	
and/or see below	

BLANK ANALYSIS RESULTS (Section 3)

Blank Actions

Qualify samples based on the criteria summarized in Table 5:

Table 5. Blank and TCLP/SPLP LEB Actions for Semivolatile Analysis

Blank Type	Blank Result	Sample Result	Action
	Detect	Non-detect	No qualification
	< CRQL	< CRQL	Report at CRQL and qualify as non-detect (U)
		> CRQI,	Use professional judgment
	**	< CRQL	Report at CRQL and qualify as non-detect (U)
Method,	> CRQL	> CRQL but < Blank Result	Report at sample results and qualify as non-detect (U) or as unusable (R)
TCLP/SPLP LEB, Field		≥ CRQL and ≥ Blank Result	Use professional judgment
	Grossly high	Detect	Report at sample results and qualify as unusable (R)
	TIC > 5.0 ug/L (water) or 0.0050 mg/L (TCLP leachate) or TIC > 170 ug/Kg (soil)	Detect	Use professional judgment

List samples qualified

CONTAMINATION SOURCE/LEVEL	COMPOUND	CONC/UNITS	AL/UNITS	SQL	AFFECTED SAMPLES

All criteria were met _X	
Criteria were not met	
and/or see below	

SURROGATE SPIKE RECOVERIES - DEUTERATED MONITORING COMPOUNDS (DMCs)

Laboratory performance of individual samples is established by evaluation of surrogate spike recoveries – deuterated monitoring compounds. All samples are spiked with surrogate compounds prior to sample analysis. The accuracy of the analysis is measured by the surrogate percent recovery. Since the effects of the sample matrix are frequently outside the control of the laboratory and may present relatively unique problems, the validation of data is frequently subjective and demands analytical experience and professional judgment.

Notes: Recoveries for DMCs in samples and blanks must be within the limits specified in Table 6.

The recovery limits for any of the compounds listed in Table 6 may be expanded at any time during the period of performance if USEPA determines that the limits are too restrictive.

If a DMC is not added in the samples and blanks or the concentrations of DMCs in the samples and blank not the specified, use professional judgment in qualifying the data.

Table 7. DMC Actions for Semivolatile Analysis

	Action		
Criteria	Detect	Non-detect	
%R < 10% (excluding DMCs with 10% as a lower acceptance limit)	.J-	R	
10% ≤ %R (excluding DMCs with 10% as a lower acceptance limit) ≤ Lower Acceptance Limit	j-	£:1	
Lower Acceptance limit < %R < Upper Acceptance Limit	No qualification	No qualification	
%R ≥ Upper Acceptance Limit	J (No qualification	

List the percent recoveries (%Rs) which do not meet the criteria for DMCs (surrogate) recovery.

Matrix:___Groundwater/Soil_____

SAMPLE ID SURROGATE COMPOUND ACTION

_DMCs_meet_the_required_criteria._Non-deuterated_surrogates_added_to_the_samples_were____
_within_laboratory_recovery_limits.______

Table 8. Semivolatile DMCs and the Associated Target Analytes

1,4-Dioxane-ds (DMC-1)	Phenol-d ₅ (DMC-2)	Bis(2-Chloroethyl) ether-d ₈ (DMC-3)
1,4-Dioxane	Benzaldehyde	Bis(2-chloroethyl)ether
	Phenol	2,2'-Oxybis(1-chloropropane)
		Bis(2-chloroethoxy)methane
2-Chlorophenol-d4 (DMC-4)	4-Methylphenol-da (DMC-5)	4-Chloroaniline-d4 (DMC-6)
2-Chlorophenol	2-Methylphenol	4-Chloroaniline
	3-Methylphenol	Hexachlorocyclopentadiene
	4-Methylphenol	Dichlorobenzidine
	2,4-Dimethylphenol	
Nitrobenzene-ds(DMC-7)	2-Nitrophenol-d4 (DMC-8)	2,4-Dichlorophenol-d3(DMC-9)
Acetophenone	1sophorone	2,4-Dichlorophenol
N-Nitroso-di-n-propylamine	2-Nitrophenol	Hexachlorobutadiene
Hexachloroethane		Hexachlorocyclopentadiene
Nitrobenzene		4-Chloro-3-methylphenol
2,6-Dinitrotoluene		2,4,6-Trichlorophenol
2,4-Dinitrotoluene		2,4.5-Trichlorophenol
N-Nitrosodiphenylamine		1,2,4,5-Tetrachlorobenzene
		*Pentachlorophenol
		2,3,4,6-Tetrachlorophenol
Dimethylphthalate-d ₆ (DMC-10)	Acenaphthylene-da(DMC-11)	4-Nitrophenol-d4(DMC-12)
Caprolactam	*Naphthalene	2-Nitroaniline
1,1'-Biphenyt	*2-Methylnaphthalene	3-Nitroaniline
Dimethylphthalate	2-Chloronaphthalene	2,4-Dinitrophenol
Diethylphthalate	*Acenaphthylene	4-Nitrophenol
Di-n-butylphthalate	*Acenaphthene	4-Nitroaniline
Butylbenzylphthalate		
Bis(2-ethylhexyl) phthalate		
Di-n-octylphthalate		1

Fluorene-d ₁₀ (DMC-13)	4,6-Dinitro-2-methylphenol-d ₂ (DMC-14)	Anthracene-d _{in} (DMC-15)	
Dibenzofuran *Fluorene 4-Chlorophenyl-phenylether 4-Bromophenyl-phenylether Carbazole	4,6-Dinitro-2-methylphenol	Hexachlorobenzene Atrazine *Phenanthrene *Anthracene	
Pyrene-d ₁₀ (DMC-16)	Benzo(a)pyrene-d ₁₂ (DMC-17)		
*l'luoranthene	3,3'-Dichlombenzidine		
*Pyrene	*Benzo(b)fluoranthene		
*Benzo(a)antliracene	*Benzo(k)fluoranthene		
*Chrysene	*Benzo(a)pyrene		
	*Indeno(1,2,3-ed)pyrene		
	*Dibenzo(a,h)anthracene		
	*Benzo(g,h,i)perylene		

^{*}Included in optional Target Analyte List (TAL) of PAHs and PCP only,

Table 9. Semivolatile SIM DMCs and the Associated Target Analytes

Fluoranthene-d10 (DMC-1)	2-Methylnaphthalene-d10 (DMC-2)		
Fluoranthene	Naphthalene		
Pyrene	2-Methylnaphthalene		
Benzo(a)anthracene	Acenaphthylene		
Chrysene	Acenaphthetic		
Benzo(b)fluoranthene	Fluorene		
Benzo(k)fluoranthene	Pentachlorophenol		
Benzo(a)pyrene	Phenanthrene		
Indeno(1,2,3-ed)pyrene	Anthracene		
Dibenzo(a,h)anthracene			
Benzo(g,h,i)perylene			

All criteria were met	
Criteria were not met	
and/or see belowX	

VII. A MATRIX SPIKE/MATRIX SPIKE DUPLICATE (MS/MSD)

This data is generated to determine long term precision and accuracy in the analytical method for various matrices. This data alone cannot be used to evaluate the precision and accuracy of individual samples. If any % R in the MS or MSD falls outside the designated range, the reviewer should determine if there are matrix effects, i.e. LCS data are within the QC limits but MS/MSD data are outside QC limit.

1. MS/MSD Recoveries and Precision Criteria

The laboratory should use one MS and a duplicate analysis of an unspiked field sample if target analytes are expected in the sample. If target analytes are not expected, MS/MSD should be analyzed.

NOTES:

Data for MS and MSDs will not be present unless requested by the Region. Notify the Contract Laboratory COR if a field or trip blank was used for the MS and MSD.

For a Matrix Spike that does not meet criteria, apply the action to only the field sample used to prepare the Matrix Spike sample. If it is clearly stated in the data validation materials that the samples were taken through incremental sampling or some other method guaranteeing the homogeneity of the sample group, then the entire sample group may be qualified.

List the %Rs, RPD of the compounds which do not meet the criteria.

Sample ID:JC21080-8A	Matrix/Level:Soil
Sample ID:JC21194-5	Matrix/Level:Soil
Sample ID:JC21223-19Q_(SIM)	Matrix/Level:Soil
The QC reported here applies to the following samples: JC21103-4, JC21103-6	Method: SW846 8270D

Compound	JC21080-8A ug/kg Q	Spike ug/kg	MS ug/kg	MS %	Spike ug/kg	MSD ug/kg	MSD %	RPD	Limits Rec/RPD
4-Bromophenyl phenyl ether 2-Chloronaphthalene	ND	9160	3850	42	9320	3200	34* c	18	39-124/33
	ND	9160	4010	44	9320	3360	36* c	18	38-110/32
4-Chlorophenyl phenyl ether 2,4-Dinitrotoluene 3,3'-Dichlorobenzidine	ND	9160	3790	41	9320	3220	35* c	16	38-119/33
	ND	9160	3360	37	9320	2530	27* c	28	28-126/36
	ND	18300	1360	7* c	18600	1490	8* c	9	10-115/44
Hexachlorobenzene Hexachlorocyclo- pentadiene Hexachloroethane 4-Nitroaniline	ND	9160	3400	37	9320	2900	31* c	16	34-125/34
	ND	18300	ND	0* c	18600	ND	0* c	nc	10-127/46
	ND	9160	1120	12* c	9320	679	7* c	49* d	21-109/38
	ND	9160	1850	20* c	9320	1650	18* c	11	21-117/38

- (a) Dilution required due to matrix interference.
- (b) Dilution required due to matrix interference. Elevated detection limit due to low volume of bad matrix sample extracted.
- (c) Outside control limits due to matrix interference.
- (d) Outside of in house control limits.
- * outside control limits

Note: MS/MSD % recoveries results and RPD apply to unspiked sample. Unspiked sample was from another project. No qualifications made.

The QC reported here applies to the following samples:

Method: SW846 8270D

JC21103-3, JC21103-5

	JC21194-5	Spike	MS	MS	Spike	MSD	MSD		Limits
Compound	ug/kg Q	ug/kg	ug/kg	%	ug/kg	ug/kg	%	RPD	Rec/RPD
Phenanthrene	1070	1950	2860	92	1990	4060	151* a	35	14-144/44

Note: MS/MSD % recoveries results and RPD apply to unspiked sample. Unspiked sample was from another project. No qualifications made.

The QC reported here applies to the following samples: JC21103-3, JC21103-4; JC21103-5, JC21103-6

Method: SW846 8270D BY SIM

	JC212	23-19Q	Spike	MS	MS	Spike	MSD	MSD		Limits
Compound	ug/kg	Q	ug/kg	ug/kg	%	ug/kg	ug/kg	%	RPD	Rec/RPD
1,4-Dioxane	ND		37.2	11.2	30* a	37.2	10.2	27* a	9	50-150/30

(a) Outside in house control limits due to matrix interference.

lote: MS/MSD % recoveries results and RPD apply to unspiked sample. Unspiked sample was from another project. No qualifications made.

No aqueous matrix MS/MSD sample analyzed in the sample batch. No action taken, blank spike/blank spike duplicate used to assess accuracy. % recoveries and RPD within laboratory control limits

- * QC limits are laboratory in-house performance criteria, LL = lower limit, UL = upper limit.
- * If QC limits are not available, use limits of 70 130 %.

Actions:

QUALITY	%R < LL	%R > UL
Positive results	J	J
Nondetects results	R	Accept

^{*} Outside control limit.

MS/MSD criteria apply only to the unspiked sample, its dilutions, and the associated MS/MSD samples:

If the % R for the affected compounds were < LL (or 70 %), qualify positive_results (J) and nondetects (UJ).

If the % R for the affected compounds were > UL (or 130 %), only qualify positive results (J). If 25 % or more of all MS/MSD %R were < LL (or 70 %) or if two or more MS/MSD %Rs were < 10%, qualify all positive results (J) and reject nondetects (R).

A separate worksheet should be used for each MS/MSD pair.

All criteria were mel _X
Criteria were not met
and/or see below

INTERNAL STANDARD PERFORMANCE

The assessment of the internal standard (IS) parameter is used to assist the data reviewer in determining the condition of the analytical instrumentation.

List the internal standard area of samples which do not meet the criteria.

DATE SAMPLE ID IS OUT IS AREA ACCEPTABLE ACTION RANGE

Internal area meets the required criteria of batch samples corresponding to this data package.

Action:

- If an internal standard area count for a sample or blank is greater than 200.0% of the area for the associated standard (opening CCV or mid-point standard from initial calibration) (see Table 10 below):
 - a. Qualify detects for compounds quantitated using that internal standard as estimated low (J-).
 - Do not qualify non-detected associated compounds.
- 2. If an internal standard area count for a sample or blank is less than 20.0% of the area for the associated standard (opening CCV or mid-point standard from initial calibration):
 - a. Qualify detects for compounds quantitated using that internal standard as estimated high (J+).
 - b. Qualify non-detected associated compounds as unusable (R).
- 3. If an internal standard area count for a sample or blank is greater than or equal to 50.0%, and less than or equal to 200% of the area for the associated standard opening CCV or mid-point standard from initial calibration, no qualification of the data is necessary.
- 4. If an internal standard RT varies by more than 10.0 seconds: Examine the chromatographic profile for that sample to determine if any false positives or negatives exist. For shifts of a large magnitude, the reviewer may consider partial or total rejection of the data for that sample fraction. Detects should not need to be qualified as unusable (R) if the mass spectral criteria are met.
- 5. If an internal standard RT varies by less than or equal to 10.0 seconds, no qualification of the data is necessary.

Note: Inform the Contract Laboratory Program Project Officer (CLP PO) if the internal standard performance criteria are grossly exceeded. Note in the Data Review Narrative potential effects on the data resulting from unacceptable internal standard performance.

State in the Data Review Narrative if the required internal standard compounds are not added to a sample or blank or if the required internal standard compound is not analyzed at the specified concentration.

Actions:

Table 10. Internal Standard Actions for Semivolatile Analysis

(Carles and a	Action			
Criteria	Detect	Non-detect		
Area response < 20% of the opening CCV or mid-point standard CS3 from ICAL	J+	R		
20% < Area response < 50% of the opening CCV or mid-point standard CS3 from ICAL	J ŧ	Ü		
50% < Area response < 200% of the opening CCV or mid-point standard CS3 from ICAL	No qualification	No qualification		
Area response > 200% of the opening CCV or mid-point standard CS3 from ICAL	J-	No qualification		
RT shift between sample/blank and opening CCV or mid-point standard CS3 from ICAL > 10.0 seconds	R	R		
RT shift between sample/blank and opening CCV or mid-point standard CS3 from ICAL < 10.0 seconds	No qualification	No qualification		

		All criteria were metX Criteria were not met and/or see below
TARGET CO	OMPOUND IDENTIFICATION	
Criteria:		
Is the Relativ RRT [openicalibration].	ve Retention Times (RRTs) of reported compoung Continuing Calibration Verification (CCV	unds within ±0.06 RRT units of the standard i) or mid-point standard from the initial Yes? or No?
List compour	nds not meeting the criteria described above:	
Sample ID	Compounds	Actions
spectrum fro	a of the sample compound and a current labor the associated calibration standard (openimust match according to the following criteria: All ions present in the standard mass spect must be present in the sample spectrum. The relative intensities of these ions must as sample spectra (e.g., for an ion with an abothe corresponding sample ion abundance more lons present at greater than 10% in the said standard spectrum, must be evaluated by interpretation.	ng CCV or mid-point standard from initial rum at a relative intensity greater than 10% gree within ±20% between the standard and undance of 50% in the standard spectrum, ust be between 30-70%). mple mass spectrum, but not present in the
List compoun	nds not meeting the criteria described above:	
Sample ID	Compounds	Actions
Identified co	ompounds meet the required criteria	

Action:

- 1. The application of qualitative criteria for GC/MS analysis of target compounds requires professional judgment. It is up to the reviewer's discretion to obtain additional information from the laboratory. If it is determined that incorrect identifications were made, qualify all such data as unusable (R).
- 2. Use professional judgment to qualify the data if it is determined that cross-contamination has occurred.
- Note in the Data Review Narrative any changes made to the reported compounds or concerns regarding target compound identifications. Note, for Contract Laboratory COR action, the necessity for numerous or significant changes.

TENTATIVELY IDENTIFIED COMPOUNDS (TICS)

NOTE: Tentatively identified compounds should only be evaluated when requested by a party from outside of the Hazardous Waste Support Section (HWSS).

Sample ID	Compound	Sample ID	Compound
	=======================================		

Action:

- 1. Qualify all TIC results for which there is presumptive evidence of a match (e.g. greater than or equal to 85% match) as tentatively identified (NJ), with approximated concentrations. TICs labeled "unknown" are qualified as estimated (J).
- 2. General actions related to the review of TIC results are as follows:
 - a. If it is determined that a tentative identification of a non-target compound is unacceptable, change the tentative identification to "unknown" or another appropriate identification, and qualify the result as estimated (J).
 - b. If all contractually-required peaks were not library searched and quantitated, the Region's designated representative may request these data from the laboratory.
- 3. In deciding whether a library search result for a TIC represents a reasonable identification, use professional judgment. If there is more than one possible match, report the result as "either compound X or compound Y". If there is a lack of isomer specificity, change the TIC result to a nonspecific isomer result (e.g., 1,3,5-trimethyl benzene to trimethyl benzene isomer) or to a compound class (e.g., 2-methyl, 3-ethyl benzene to a substituted aromatic compound).
- 4. The reviewer may elect to report all similar compounds as a total (e.g., all alkanes may be summarized and reported as total hydrocarbons).

- 5. Target compounds from other fractions and suspected laboratory contaminants should be marked as "non-reportable".
- 6. Other Case factors may influence TIC judgments. If a sample TIC match is poor, but other samples have a TIC with a valid library match, similar RRT, and the same ions, infer identification information from the other sample TIC results.
- 7. Note in the Data Review Narrative any changes made to the reported data or any concerns regarding TIC identifications.
- 8. Note, for Contract Laboratory COR action, failure to properly evaluate and report TICs

All criteria were met _	_X
Criteria were not met	100
and/or see below	

SAMPLE QUANTITATION AND REPORTED CONTRACT REQUIRED QUANTITATION LIMITS (CRQLS)

Action:

- 1. When a sample is analyzed at more than one dilution, the lower CRQL are used unless a QC exceedance dictates the use of higher CRQLs from the diluted sample. Samples reported with an "E" qualifier should be reported from the diluted sample.
- 2. If any discrepancies are found, the Region's designated representative may contact the laboratory to obtain additional information that could resolve any differences. If a discrepancy remains unresolved, the reviewer must use professional judgment to decide which value is the most accurate. Under these circumstances, the reviewer may determine that qualification of data is warranted. Note in the Data Review Narrative a description of the reasons for data qualification and the qualification that is applied to the data.
- 3. For non-aqueous samples, if the solids is less than 10.0%, use professional judgment for both detects and non-detects. If the percent solid for a soil sample is greater than or equal to 10.0% and less than 30.0%, use professional judgment to qualify detects and non-detects. If the percent solid for a soil sample is greater than or equal to 30.0%, detects and non-detects should not be qualified (see Table 11).
- 4. Note, for Contract Laboratory COR action, numerous or significant failures to accurately quantify the target compounds or to properly evaluate and adjust CRQLs.
- 5. Results between MDL and CRQL should be qualified as estimated "J".
- 6. Results < MDL should be reported at the CRQL and qualified "U". MDLs themselves should not be reported.

Table 11. Percent Solids Actions for Semivolatile Analysis for Non-Aqueous Samples

(Suite main	Ac	Action				
Criteria	Detects	Non-detects				
%Solids < 10.0%	Use professional judgment	Use professional judgment				
10.0% < %Solids < 30.0%	Use professional judgment	Use professional judgment				
%Solids > 30.0%	No qualification ,	No qualification				

SAMPLE QUANTITATION

The sample quantitation evaluation is to verify laboratory quantitation results. In the space below, please show a minimum of one sample calculation:

QUANTITATION LIMITS

A. Dilution performed

SAMPLE ID	DILUTION FACTOR	REASON FOR DILUTION	
	· ·		 _

	All criteria were metN/A Criteria were not met and/or see below
FIELD DUPLICATE PRECISION	
Sample IDs:	Matrix:

Field duplicates samples may be taken and analyzed as an indication of overall precision. These analyses measure both field and lab precision; therefore, the results may have more variability than laboratory duplicates which only laboratory performance. It is also expected that soil duplicate results will have a greater variance than water matrices due to difficulties associated with collecting identical field duplicate samples.

The project QAPP should be reviewed for project-specific information.

Suggested criteria: if large RPD (> 50 %) is observed, confirm identification of the samples and note differences. If both samples and duplicate are <5 SQL, the RPD criteria is doubled.

COMPOUND	SQL ug/L	SAMPLE CONC.	DUPLICATE CONC.	RPD	ACTION
No field/laborator used to assess pr	ry duplicate recision. R	e analyzed as p PD within the re	art of this data packa	nge. MS/MS for detected	D % recoveries RPD target analytes.

All criteria were metX Criteria were not met and/or see below

OTHER ISSUES

A. Sys	tem Perfo	ormance
--------	-----------	---------

List samples qualified based on	the degradation of system	performance during simple analysi
---------------------------------	---------------------------	-----------------------------------

Sample ID	Comments	Actions			
	-				

Action:

Use professional judgment to qualify the data if it is determined that system performance has degraded during sample analyses. Inform the Contract Laboratory Program COR any action as a result of degradation of system performance which significantly affected the data.

B. Overall Assessment of Data

List samples qualified based on other issues:

Sample ID	Comments	Actions			
No_other_issues_ for_decission_pur		_the_dataResults_are_valid_and_can_be_use	ed		

Action:

- 1. Use professional judgment to determine if there is any need to qualify data which were not qualified based on the Quality Control (QC) criteria previously discussed.
- 2. Write a brief narrative to give the user an indication of the analytical limitations of the data. Inform the Contract Laboratory COR the action, any inconsistency of the data with the Sample Delivery Group (SDG) Narrative. If sufficient information on the intended use and required quality of the data is available, the reviewer should include their assessment of the usability of the data within the given context. This may be used as part of a formal Data Quality Assessment (DQA).
- 3. Sometimes, due to dilutions, re-analysis or SIM/Scan runs are being performed, there will be multiple results for a single analyte from a single sample. The following criteria and professional judgment are used to determine which result should be reported:
 - The analysis with the lower CRQL
 - The analysis with the better QC results
 - The analysis with the higher results

EXECUTIVE NARRATIVE

SDG No:

JC21103

Laboratory:

Accutest, New Jersey

Analysis:

SW846-8081B

Number of Samples:

2

Location:

BMSMC, Building 5 Area

Humacao, PR

SUMMARY:

Two (2) sample was analyzed for selected pesticides following method SW846-8081B. The sample results were assessed according to USEPA data validation guidance documents in the following order of precedence *Hazardous Waste Support Section SOP No. HW-36A, Revision O, June, 2015. SOM02.2. Pesticide Data Validation.* The QC criteria and data validation actions listed on the data review worksheets are from the primary suidance document upless otherwise peted.

guidance document, unless otherwise noted.

Results are valid and can be used for decision making purposes.

Critical issues:

None

Major:

None

Minor:

None

Critical findings:

None

Major findings:

None

Minor findings:

1. No MS/MSD analyzed with this data package. Blank spike/blank spike % recoveries used to assess accuracy. % recoveries and RPD within laboratory control limits. No action

taken.

COMMENTS:

Results are valid and can be used for decision making purposes.

Reviewers Name:

Rafael Infante

Chemist License 1888

Signature:

Date:

June 21, 2016

SAMPLE ORGANIC DATA SAMPLE SUMMARY

Sample ID: JC21103-1

Sample location: BMSMC Building 5 Area

Sampling date: 25-May-16 Matrix: Groundwater

METHOD: 8081B

METHOE), OUOID					
Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Aldrin	0.011	ug/L	1	-	U	Yes
alpha-BHC	0.011	ug/L	1	-	U	Yes
beta-BHC	0.011	ug/L	1	-	U	Yes
delta-BHC	0.011	ug/L	1	-	U	Yes
gamma-BHC (Lindane)	0.011	ug/L	1	-	U	Yes
alpha-Chlordane	0.011	ug/L	1	-	U	Yes
gamma-Chlordane	0.011	ug/L	1	-	U	Yes
Dieldrin	0.011	ug/L	1	-	U	Yes
4,4'-DDD	0.011	ug/L	1	-	บ	Yes
4,4'-DDE	0.011	ug/L	1	•	U	Yes
4,4'-DDT	0.011	ug/L	1		U	Yes
Endrin	0.011	ug/L	1	-	U	Yes
Endosulfan sulfate	0.011	ug/L	1	-	U	Yes
Endrin aldehyde	0.011	ug/L	1	-	U	Yes
Endrin ketone	0.011	ug/L	1	-	U	Yes
Endosulfan-I	0.011	ug/L	1	12	U	Yes
Endosulfan-II	0.011	ug/L	1	22	U	Yes
Heptachlor	0.011	ug/L	1	2	U	Yes
Heptachlor epoxide	0.011	ug/L	1	1.7	U	Yes
Methoxychlor	0.022	ug/L	1	12	y U	Yes
Toxaphene	0.28	ug/L	1	17	U	Yes

Sample ID: JC21103-2

Sample location: BMSMC Building 5 Area

Sampling date: 25-May-16

Matrix: AQ - Equipment Blank

METHOD: 8081B

Analyte Name	Re	sult	Units	Dilution Factor	Lab Flag	Validation	Reportable
Aldrin	0.0)10	ug/L	1	-	U	Yes
alpha-BHC	0.0)10	ug/L	1	-	U	Yes
beta-BHC	0.0)10	ug/L	1	-	U	Yes
delta-BHC	0.0	010	ug/L	1	-	U	Yes
gamma-BHC (Lindane)	0.0	010	ug/L	1	-	U	Yes
alpha-Chlordane	0.0)10	ug/L	1	-	U	Yes
gamma-Chlordane	0.0)10	ug/L	1	-	U	Yes
Dieldrin	0.0)10	ug/L	1	-	U	Yes
4,4'-DDD	0.0)10	ug/L	1	-	U	Yes
4,4'-DDE	0.0)10	ug/L	1	-	Ü	Yes
4,4'-DDT	0.0)10	ug/L	1	-	U	Yes
Endrin	0.0)10	ug/L	1	-	U	Yes
Endosulfan sulfate	0.0)10	ug/L	1	-	U	Yes
Endrin aldehyde	0.0)10	ug/L	1	-	U 🥫	Yes
Endrin ketone	0.0)10	ug/L	1	-	U	Yes
Endosulfan-I	0.0)10	ug/L	1	-	U	Yes
Endosulfan-II	0.0)10	ug/L	1	-	U	Yes
Heptachlor	0.0)10	ug/L	1	-	IJ	Yes
Heptachlor epoxide	0.0)10	_ ug/L	1	-	U	Yes
Methoxychlor	0.0)20	ug/L	1	-	U	Yes
Toxaphene	0.	25	ug/L	1	-	U	Yes

	Sampling Date:_	May_25-26,_2016_
	Shipping Date:	May_26,_2016
		2
REVIEW OF PESTICIDE ORG	GANIC PACKAGE	
The following guidelines for evaluating volatile required validation actions. This document will as udgment to make more informed decision and it users. The sample results were assessed according to the sample results are supplied to the sample	ssist the reviewer n better serving ting to USEPA datalazardous Waste Sude Data Validation. ew worksheets a	in using professionathe needs of the data a validation guidance apport Section SOP No The QC criteria and are from the primary
The hardcopied (laboratory name) _Accutesteviewed and the quality control and performance data summer	data padarized. The data revi	ckage received has been ew for VOCs included:
ab. Project/SDG No.:JC21103 No. of Samples:2	Sample matrix:	Groundwater
Field blank No.:		
Equipment blank No.:MC21103-2		
Field duplicate No.:		
Field enikes No:		
Field spikes No.:		
20 audit samples		<u>.</u>
X Data CompletenessX Holding TimesN/A GC/MS TuningX Internal Standard PerformanceX BlanksX Surrogate RecoveriesX Matrix Spike/Matrix Spike Duplicate Overall Comments:TCL_pesticides_list_by_SW846-8	X Field DuX CalibratX CompouX CompouX Quantita	ions und Identifications und Quantitation
Definition of Qualifiers: - Estimated results U- Compound not detected R- Rejected data UJ- Estimated nandetect Reviewer: Au Man		
Date:June_22,_2016		

Project/Case Number:_____JC21103____

DATA COMPLETENESS

MISSING INFORMATION	DATE LAB. CONTACTED	DATE RECEIVED
		_ ···
		<u> </u>
<u> </u>	2000	
		WWW 7—17 —

All criteria were metX
Criteria were not met
and/or see below

HOLDING TIMES

The objective of this parameter is to ascertain the validity of the results based on the holding time of the sample from time of collection to the time of analysis.

Complete table for all samples and note the analysis and/or preservation not within criteria

SAMPLE ID	DATE SAMPLED	DATE EXTRACTED/ANALYZED	ACTION
Samples property pro	eserved.		
	<u> </u>		

Preservatives:	_All_samples_extrac	ted_and_analyzed_wi	thin_the_required_	criteria

<u>Criteria</u>

Aqueous samples - seven (7) days from sample collection for extraction; 40 days from sample collection for analysis.

Non-aqueous samples – fourteen (14) days from sample collection for extraction; 40 days from sample collection for analysis.

Cooler temperature (Criteria: 4 ± 2 °C): 4°C - OK

Actions

Qualify aqueous sample results using preservation and technical holding time information as follows:

- a. If there is no evidence that the samples were properly preserved ($T = 4^{\circ}C \pm 2^{\circ}C$), and the samples were extracted or analyzed within the technical holding times, qualify detects as estimated (J) and non-detects as estimated (UJ).
- b. If there is no evidence that the samples were properly preserved (T = 4° C \pm 2° C), and the samples were extracted or analyzed outside the technical holding times, qualify detects as estimated (J) and non-detects as estimated (UJ).
- c. If the samples were properly preserved, and were extracted and analyzed within the technical holding times, no qualification of the data is necessary.
- d. If the samples were properly preserved, and were extracted or analyzed outside the technical holding times, qualify detects as estimated (J) and non-detects as estimated (UJ). Note in the Data Review Narrative that holding times were exceeded and the effect of exceeding the holding time on the resulting data.

- e. Use professional judgment to qualify samples whose temperature upon receipt at the laboratory is either below 2 degrees centigrade or above 6 degrees centigrade.
- f. If technical holding times are grossly exceeded, use professional judgment to qualify the data.

Qualify non-aqueous sample results using preservation and technical holding time information as follows:

- a. If there is no evidence that the samples were properly preserved (T = 4° C \pm 2° C), and the samples were extracted or analyzed within the technical holding time, qualify detects as estimated (J) and non-detects as estimated (UJ).
- b. If there is no evidence that the samples were properly preserved (T = 4° C \pm 2° C), and the samples were extracted or analyzed outside the technical holding time, qualify detects as estimated (UJ).
- c. If the samples were properly preserved, and were extracted and analyzed within the technical holding time, no qualification of the data is necessary.
- d. If the samples were properly preserved, and were extracted or analyzed outside the technical holding time, qualify detects as estimated (J) and non-detects as estimated (UJ). Note in the Data Review Narrative that holding times were exceeded and the effect of exceeding the holding time on the resulting data.
- e. Use professional judgment to qualify samples whose temperature upon receipt at the laboratory is either below 2 degrees centigrade or above 6 degrees centigrade.
- f. If technical holding times are grossly exceeded, use professional judgment to qualify the data.

All criteria were met _	_x_	_
Criteria were not met see belo	W	

GAS CHROMATOGRAPH WITH ELECTRON CAPTURE DETECTOR (GC/ECD) INSTRUMENT PERFORMANCE CHECK (SECTIONS 1 TO 5)

1. Resolution Check Mixture

Criteria

Is the resolution between two adjacent peaks in the Resolution Check Mixture C greater than or equal to 80.0% for all analytes for the primary column and greater than or equal to 50.0% for the confirmation column?

Yes? or No?

Is the resolution between two adjacent peaks in the Resolution Check Mixture (A and B) greater than or equal to 60.0%?

Yes? or No?

Note:

If resolution criteria are not met, the quantitative results may not be accurate due to inadequate resolution. Qualitative identifications may also be questionable if coelution exists.

Action

- a. Qualify detects for target compounds that were not adequately resolved as tentatively identified (NJ).
- b. Qualify non-detected compounds as unusable (R).

2. Performance Evaluation Mixture (PEM) Resolution Criteria

Criteria

Is PEM analysis performed at the required frequency (at the end of each pesticide initial calibration sequence and every 12 hours)?

Yes? or No?

Action

a. If PEM is not performed at the required frequency, qualify all associated sample and blank results as unusable (R).

Criteria

Is PEM % Resolution < 90%?

Yes? or No?

Action

- a. a. Qualify detects for target compounds that were not adequately resolved as tentatively identified (NJ).
- b. Qualify non-detected compounds as unusable (R).

All criteria were met ___X__
Criteria were not met see below ____

3. PEM 4,4'-DDT Breakdown

Criteria

Is the PEM 4,4'-DDT % Breakdown >20.0% and 4,4'-DDT is detected?

Yes? or No?

Action

a. Qualify detects for 4,4'-DDT; detects for 4,4'-DDD; and detects for 4,4'-DDE as estimated (J)

Criteria

is the PEM 4,4'-DDT % Breakdown >20.0% and 4,4'-DDT is not detected

Yes? or No?

Action

- a. Qualify non-detects for 4,4'- DDT as unusable (R)
- b. Qualify detects for 4,4'-DDD as tentatively identified (NJ)
- c. Qualify detects for 4,4'-DDE as tentatively identified (NJ)

4. PEM Endrin Breakdown

Criteria

Is the PEM Endrin % Breakdown >20.0% and Endrin is detected?

Yes? or No?

Action

a. Qualify detects for Endrin; detects for Endrin aldehyde; and detects for Endrin ketone as estimated (J)

Criteria

Is the PEM Endrin % Breakdown >20.0% and Endrin is not detected

Yes? or No?

Action

- a. Qualify non-detects for Endrin as unusable (R)
- b. Qualify detects for Endrin aldehyde as tentatively identified (NJ)
- c. Qualify detects for Endrin ketone as tentatively identified (NJ)

All criteria were met _	_X
Criteria were not met see belo	W

5. Mid-point Individual Standard Mixture Resolution -

Criteria

Is the resolution between two adjacent peaks in the Resolution Check Mixture C greater than or equal to 80.0% for all analytes for the primary column and greater than or equal to 50.0% for the confirmation column?

Yes? or No?

Is the resolution between two adjacent peaks in the Resolution Check Mixture (A and B) greater than or equal to 90.0%?

Yes? or No?

Note: If resolution criteria are not met, the quantitative results may not be accurate due to inadequate resolution. Qualitative identifications may also be questionable if coelution exists.

Action

- a. Qualify detects for target compounds that were not adequately resolved as tentatively identified (NJ).
- b. Qualify non-detected compounds as unusable (R).

Criteria

Is mid-point individual standard mixture analysis performed at the required frequency (every 12 hours)?

Yes? or No?

Action

a. If the mid-point individual standard mixture analysis is not performed at the required frequency, qualify all associated sample and blank results as unusable (R).

All criteria were metX Criteria were not met and/or see below
stablished to ensure that the ative data.
05/24/16
fication:_05/24/16
n:_06/01/16
06/02/16
GC4G
Aqueous/low
•
SAMPLES AFFECTED
-
 ent performance criteria. in at least one of the two leet the performance criteria
s shown in Table 3 of SOF <u>Yes</u> ? or No?
onal judgment to evaluate the

CALIBRATION VERIFICATION

Compliance requirements for satisfactory instrument calibration are established to ensure that the instrument is capable of producing and maintaining acceptable quantitative data.

Date of initial calibration:	05/24/16
Dates of initial calibration v	rerification:_05/24/16
Dates of continuing calibra	tion:_06/01/16
Dates of final calibration	06/02/16
Instrument ID numbers:	GC4G
Matrix/Level:	Aqueous/low

DATE	LAB	FILE	CRITERIA OUT	COMPOUND	SAMPLES AFFECTED
	ID#		RFs, %RSD, %D, r		
	•			<u> </u>	
,	T				
Initia	l and init	ial calib	ration verification within	the guidance docu	ment performance criteria.
					ria in at least one of the two
columns.	Final ca	libration	verification included in	data. % differences	s meet the performance criteria
				the two columns.	
		_			
	+			-	
			L	<u> </u>	

Criteria

Are a five point calibration curve delivered with concentration levels as shown in Table 3 of SOP HW-36A, Revision 0, June, 2015?

Yes? or No?

Actions

If the standard concentrations listed in Table 3 are not used, use professional judgment to evaluate the effect on the data

Criteria

Are RT Windows calculated correctly?

Yes? or No?

Action

Recalculate the windows and use the corrected values for all evaluations.

Criteria

Are the Percent Relative Standard Deviation (%RSD) of the CFs for each of the single component target compounds less than or equal to 20.0%, except for alpha-BHC and delta-BHC?

Yes? or No?

Are the %RSD of the CFs for alpha-BHC and delta-BHC less than or equal to 25.0%. Yes? or No?

Is the %RSD of the CFs for each of the Toxaphene peaks must be < 30% when 5-point ICAL is performed?

Yes? or No?

Is the %RSD of the CFs for the two surrogates (tetrachloro-m-xylene and decachlorobiphenyl) less than or equal to 30.0%.

Yes? or No?

Action

- a. If the %RSD criteria are not met, qualify detects as estimated (J) and use professional judgment to qualify non-detected target compounds.
- b. If the %RSD criteria are within allowable limits, no qualification of the data is necessary

Continuing Calibration Checks

Criteria

Is the continuing calibration standard analyzed at the acceptable time intervals? Yes? or No?

Action

- a. If more than 14 hours has elapsed from the injection of the instrument blank that begins an analytical sequence (opening CCV) and the injection of either a PEM or mid-point concentration of the Individual Standard Mixtures (A and B) or (C), qualify all data as unusable (R).
- b. If more than 12 hours has elapsed from the injection of the instrument blank that begins an analytical sequence (opening CCV) and the injection of the last sample or blank that is part of the same analytical sequence, qualify all data as unusable (R).
- c. If more than 72 hours has elapsed from the injection of the sample with a Toxaphene detection and the Toxaphene Calibration Verification Standard (CS3), qualify all data as unusable (R).

Criteria

Is the Percent Difference (%D) within ±25.0% for the PEM sample?

Yes? or No?

Action

a. Qualify associated detects as estimated (J) and non-detects as estimated (UJ).

Criteria

For the Calibration Verification Standard (CS3); is the Percent Difference (%D) within ±25.0%? Yes? or No?

Action

Qualify associated detects as estimated (J) and non-detects as estimated (UJ).

All criteria were met _	х_
Criteria were not met	
and/or see below	

Criteria

Is the PEM 4,4'-DDT % Breakdown >20.0% and 4,4'-DDT is detected?

Yes? or No?

Action

- a. Qualify detects for 4,4'-DDT; detects for 4,4'-DDD; and detects for 4,4'-DDE as estimated (J)
- b. Non-detected associated compounds are not qualified

Criteria

Is the PEM 4,4'-DDT % Breakdown >20.0% and 4,4'-DDT is not detected

Yes? or No?

Action

- a. Qualify non-detects for 4,4'- DDT as unusable (R)
- b. Qualify detects for 4,4'-DDD as tentatively identified (NJ)
- c. Qualify detects for 4,4'-DDE as tentatively identified (NJ)

Criteria

Is the PEM Endrin % Breakdown >20.0% and Endrin is detected?

Yes? or No?

Action

- a. Qualify detects for Endrin; detects for Endrin aldehyde; and detects for Endrin ketone as estimated (J)
- b. Non-detected associated compounds are not qualified

Criteria

Is the PEM Endrin % Breakdown >20.0% and Endrin is not detected

Yes? or No?

Action

- a. Qualify non-detects for Endrin as unusable (R)
- b. Qualify detects for Endrin aldehyde as tentatively identified (NJ)
- c. Qualify detects for Endrin ketone as tentatively identified (NJ)

A separate worksheet should be filled for each initial curve

All criteria were metX
Criteria were not met
and/or see below

BLANK ANALYSIS RESULTS (Sections 1 & 2)

The assessment of the blank analysis results is to determine the existence and magnitude of contamination problems. The criteria for evaluation of blanks apply only to blanks associated with the samples, including trip, equipment, and laboratory blanks. If problems with any blanks exist, all data associated with the case must be carefully evaluated to determine whether or not there is an inherent variability in the data for the case, or if the problem is an isolated occurrence not affecting other data.

List the contami	ination in the bla	anks below. Hig	th and low levels blanks	must be treated separately.
CRQL concentr	ationN	/A		
Laboratory blan	ks			
DATE Analyzed	LABID	LEVEL/ MATRIX	COMPOUND	CONCENTRATION UNITS
			왕님 [19] [1]	it_of_0.01_and_0.001_ug/L
Field/ <u>Equipmen</u>				
	LABID	LEVEL! MATRIX	COMPOUND	CONCENTRATION UNITS
	alytes_detected_	MATRIX in_the_equipm		UNITS blanks_analyzed_with_this
ANALYZED _No_target_ana _data_package.	alytes_detected_	MATRIX _in_the_equipm	nent_blankNo_field/trip	UNITS _blanks_analyzed_with_this
ANALYZED _No_target_ana _data_package.	alytes_detected_	MATRIX _in_the_equipn	nent_blankNo_field/trip	UNITS _blanks_analyzed_with_this
ANALYZED _No_target_ana _data_package.	alytes_detected_	MATRIX _in_the_equipn	nent_blankNo_field/trip	UNITS _blanks_analyzed_with_this

All criteria were met	_X_	
Criteria were not met		
and/or see below	_	

BLANK ANALYSIS RESULTS (Section 3)

Blank Actions

Action Levels (ALs) should be based upon the highest concentration of contaminant determined in any blank. Do not qualify any blank with another blank. The ALs for samples which have been diluted should be corrected for the sample dilution factor and/or % moisture, where applicable. No positive sample results should be reported unless the concentration of the compound in the samples exceeds the ALs:

The concentration of non-target compounds in all blanks must be less than or equal to 10 μ g/L. The concentration of each target compound found in the method or field blanks must be less than its CRQL listed in the method.

Data concerning the field blanks are not evaluated as part of the CCS process. If field blanks are present, the data reviewer should evaluate this data in a similar fashion as the method blanks.

Specific actions are as follows:

Blank Actions for Pesticide Analyses

Blank Type	Blank Result	Sample Result	Action for Samples
	Detects	Not detected	No qualification required
Method, Sulfur Cleanup, Instrument, Field, TCLP/SPLP	< CRQL	< CRQL	Report CRQL value with a U
		≥CRQL	No qualification required
		< CRQL	Report CRQL value with a U
	> CRQL	≥ CRQL and ≤ blank concentration	Report blank value for sample concentration with a U
		≥ CRQL and > blank concentration	No qualification required
	= CRQL	≤CRQL	Report CRQL value with a U
		> CRQL	No qualification required
	Gross contamination	Detects	Report blank value for sample concentration with a U

All criteria were metX
Criteria were not met
and/or see below

CONTAMINATION SOURCE/LEVEL	COMPOUND	CONC/UNITS	AL/UNITS	SQL	AFFECTED SAMPLES

13

All criteria were melX
Criteria were not met
and/or see below

SURROGATE SPIKE RECOVERIES

Laboratory performance of individual samples is established by evaluation of surrogate spike recoveries. All samples are spiked with surrogate compounds prior to sample analysis. The accuracy of the analysis is measured by the surrogate percent recovery. Since the effects of the sample matrix are frequently outside the control of the laboratory and may present relatively unique problems, the validation of data is frequently subjective and demands analytical experience and professional judgment.

List the percent recoveries (%Rs) which do not meet the criteria for surrogate recovery.

Matrix:_Aqueous						
Lab Sample ID	Lab File ID	S1 a	S1 b	S2 a	S2 b	
JC21103-1	4G68800.D	76	84	40	37	
JC21103-2	4G68801.D	62	59	29	27	
OP94371-BS1	4G68796.D	97	88	97	89	
OP94371-BSD	4G68797.D	88	83	79	76	
OP94371-MB1	4G68795.D	95	89	78	72	
Surrogate		Recovery				
Compounds		Limits				
S1 = Tetrachloro-m-xylene		26-132%				
S2 = Decachlorobiphenyl		10-118%				
(a) Recovery from GC signal #1 (b) Recovery from GC signal #2						

Note: Surrogate recoveries within laboratory control limits.

Actions:

- a. For any surrogate recovery greater than 150%, qualify detected target compounds as biased high (J+).
- b. Do not qualify non-detected target compounds for surrogate recovery > 150 %.
- c. If both surrogate recoveries are greater than or equal to 30% and less than or equal to 150%, no qualification of the data is necessary.
- d. For any surrogate recovery greater than or equal to 10% and less than 30%, qualify detected target compounds as biased low (J-).
- e. For any surrogate recovery greater than or equal to 10% and less than 30%, qualify non-detected target compounds as approximated (UJ).
- f. If low surrogate recoveries are from sample dilution, professional judgment should be used to determine if the resulting data should be qualified. If sample dilution is not a factor:
 - i. Qualify detected target compounds as biased low (J-).
 - ii. Qualify non-detected target compounds as unusable (R).

- g. If surrogate RTs in PEMs, Individual Standard Mixtures, samples, and blanks are outside of the RT Windows, the reviewer must use professional judgment to qualify data.
- h. If surrogate RTs are within RT windows, no qualification of the data is necessary.
- i. If the two surrogates were not added to all samples, MS/MSDs, standards, LCSs, and blanks, use professional judgment in qualifying data as missing surrogate analyte may not directly apply to target analytes.

Summary Surrogate Actions for Pesticide Analyses

	Action*			
Criteria	Detected Target	Non-detected Target		
	Compounds	Compounds		
%R > 150%	J+	No qualification		
30% < %R < 150%	No qualification			
10% < %R < 30%	J-	UJ		
%R < 10% (sample dilution not a factor)	J-	R		
%R < 10% (sample dilution is a factor)	Use professional judgment			
RT out of RT window	Use professional judgment			
RT within RT window	No qualification			

Use professional judgment in qualifying data, as surrogate recovery problems may not directly apply to target analytes.

All criteria were met	_N/A
Criteria were not met	
and/or see below	

MATRIX SPIKE/MATRIX SPIKE DUPLICATE (MS/MSD)

This data is generated to determine long term precision and accuracy in the analytical method for various matrices. This data alone cannot be used to evaluate the precision and accuracy of individual samples. If any % R in the MS or MSD falls outside the designated range, the reviewer should determine if there are matrix effects, i.e. LCS data are within the QC limits but MS/MSD data are outside QC limit.

MS/MSD Recoveries and Precision Criteria

List the %Rs. RPD of the compounds which do not meet the criteria

Data for MS and MSDs will not be present unless requested by the Region.

Notify the Contract Laboratory Program Project Officer (CLP PO) if a field blank was used for the MS and MSD, unless designated as such by the Region.

NOTE: For a Matrix Spike that does not meet criteria, apply the action to only the field sample used to prepare the Matrix Spike sample. If it is clearly stated in the data validation materials that the samples were taken through incremental sampling or some other method guaranteeing the homogeneity of the sample group, then the entire sample group may be qualified.

Sample ID:			Matrix/	/Level:	
MS OR MSD	COMPOUND	% R	RPD	QC LIMITS	ACTION
			- 28		
	3 2 2 3 3 3 3 3 3 3 3 3 3 3 3 3 3 3 3 3	150.			

Note: No MS/MSD sample analyzed with this data package. Blank spike/blank spike duplicate used to assess accuracy. % recoveries and RPD within laboratory control limits. No action taken.

Action

No qualification of the data is necessary on MS and MSD data alone. However, using professional judgment, the validator may use the MS and MSD results in conjunction with other QC criteria and determine the need for some qualification of the data.

A separate worksheet should be used for each MS/MSD pair.

All criteria were met _X
Criteria were not met
and/or see below

LABORATORY CONTROL SAMPLE (LCS) ANALYSIS

This data is generated to determine accuracy of the analytical method for various matrices.

1. LCS Recoveries Criteria

LCS Spike Compound	Recovery Limits (%)
gamma-BHC	50 – 120
Heptachlor epoxide	50 – 150
Dieldrin	30 – 130
4,4'-DDE	50 – 150
Endrin	50 – 120
Endosulfan sulfate	50 – 120
trans-Chlordane	30 – 130
Tetrachioro-m-xylene (surrogate)	30 – 150
Decachlorobiphenyl (surrogate)	30 – 150

LC3	S concentrations	:0.25_ug/l		
List the %R	of compounds v	which do not meet the criteria	ì	
	LCS ID	COMPOUND	% R	QC LIMIT
				

Action

The following guidance is suggested for qualifying sample data for which the associated LCS does not meet the required criteria.

- a. If the LCS recovery exceeds the upper acceptance limit, qualify detected target compounds as estimated (J). Do not qualify non-detected target compounds.
- b. If the LCS recovery is less than the lower acceptance limit, qualify detected target compounds as estimated (J) and non-detects as unusable (R).
- c. Use professional judgment to qualify data for compounds other than those compounds that are included in the LCS.
- d. Use professional judgment to qualify non-LCS compounds. Take into account the compound class, compound recovery efficiency, analytical problems associated with each compound, and comparability in the performance of the LCS compound to the non-LCS compound.
- e. If the LCS recovery is within allowable limits, no qualification of the data is necessary.

2. Frequency Criteria:

Where LCS analyzed at the required frequency and for each matrix? <u>Yes</u> or No. If no, the data may be affected. Use professional judgment to determine the severity of the effect and qualify data accordingly. Discuss any actions below and list the samples affected.

Note: Blank spike/blank spike duplicate analyzed for aqueous matrices. % recoveries and RPD within laboratory control limits.

All criteria were met	
Criteria were not met	
and/or see belowN/A	_

FLORISIL CARTRIDGE PERFORMANCE CHECK

NOTE: Florisil cartridge cleanup is mandatory for all extracts.

Criteria

Is the Florisil cartridge performance check conducted at least once on each lot of cartridges used for sample cleanup or every 6 months, whichever is most frequent?

Yes? or No?

Criteria

Are the results for the Florisil Cartridge Performance Check solution included with the data package?

Yes? or No?

Note: If % criteria are not met, examine the raw data for the presence of polar interferences and use professional judgment in qualifying the data as follows:

Action:

- a. If the Percent Recovery is greater than 120% for any of the pesticide target compounds in the Florisil Cartridge Performance Check, qualify detected compounds as estimated (J). Do not qualify non-detected target compounds.
- b. If the Percent Recovery is greater than or equal to 80% and less than or equal to 120% for all the pesticide target compounds, no qualification of the data is necessary.
- c. If the Percent Recovery is greater than or equal to 10% and less than 80% for any of the pesticide target compounds in the Florisil Cartridge Performance Check, qualify detected target compounds as estimated (J) and non-detected target compounds as approximated (UJ).
- d. If the Percent Recovery is less than 10% for any of the pesticide target compounds in the Florisil Cartridge Performance Check, qualify detected compounds as estimated (J) and qualify non-detected target compounds as unusable (R).
- e. If the Percent Recovery of 2,4,5-trichlorophenol in the Florisil Cartridge Performance Check is greater than or equal to 5%, use professional judgment to qualify detected and non-detected target compounds, considering interference on the sample chromatogram.

Note: State in the Data Review Narrative potential effects on the sample data resulting from the Florisil Cartridge Performance Check analysis not yielding acceptable results.

Note:_ No information for florisil cartridge performance check included in data package. Florisil cartridge was not used for sample extraction/clean-up. No qualification of the data performed, professional judgment.

All criteria were metN/A	
Criteria were not met	
and/or see below	

GEL PERMEATION CHROMATOGRAPHY (GPC) PERFORMANCE CHECK

NOTE: GPC cleanup is mandatory for all soil samples.

If GPC criteria are not met, examine the raw data for the presence of high molecular weight contaminants; examine subsequent sample data for unusual peaks; and use professional judgment in qualifying the data. Notify the Contract Laboratory Program Project Officer (CLP PO) if the laboratory chooses to analyze samples under unacceptable GPC criteria.

Action:

- a. If the Percent Recovery is less than 10% for the pesticide compounds and surrogates during the GPC calibration check, the non-detected target compounds may be suspect, qualify detected compounds as estimated (J).
- b. If the Percent Recovery is less than 10% for the pesticide compounds and surrogates during the GPC calibration check, qualify all non-detected target compounds as unusable (R).
- c. If the Percent Recovery is greater than or equal to 10% and is less than 80% for any of the pesticide target compounds in the GPC calibration, qualify detected target compounds as estimated (J) and non-detected target compounds as approximated (UJ).
- d. If the Percent Recovery is greater than or equal to 80% and less than or equal to 120% for all the pesticide target compounds, no qualification of the data is necessary.
- e. If high recoveries (i.e., greater than 120%) were obtained for the pesticides and surrogates during the GPC calibration check, qualify detected compounds as estimated (J). Do not qualify non-detected target compounds.

Note: State in the Data Review Narrative potential effects on the sample data resulting from the GPC cleanup analyses not yielding acceptable results.

Note: No information for performance of GPC cleanup included in data package. No qualification of the data performed, professional judgment.

All criteria were met	_X
Criteria were not met	
and/or see below	_

TARGET COMPOUND IDENTIFICATION

Criteria:

- 1. Is Retention Times (RTs) of both of the surrogates and reported target compounds in each sample within the calculated RT Windows on both columns?

 Yes? or No?
- 2. Is the Tetrachloro-m-xylene (TCX) RT ±0.05 minutes of the Mean RT (RT) determined from the initial calibration and Decachlorobiphenyl (DCB) within ±0.10 minutes of the RT determined from the initial calibration?

 Yes? or No?
- 3. Is the Percent Difference (%D) for the detected mean concentrations of a pesticide target compound between the two Gas Chromatograph (GC) columns within the inclusive range of \pm 25.0 %?

 Yes? or No?
- 4. When no analytes are identified in a sample; are the chromatograms from the analyses of the sample extract and the low-point standard of the initial calibration associated with those analyses on the same scaling factor?

 Yes? or No?
- 5. Does the chromatograms display the Single Component Pesticides (SCPs) detected in the sample and the largest peak of any multi-component analyte detected in the sample at less than full scale.

 Yes? or No?
- 6. If an extract is diluted; does the chromatogram display SCPs peaks between 10-100% of full scale, and multi-component analytes between 25-100% of full scale? Yes? or No?
- 7. For any sample; does the baseline of the chromatogram return to below 50% of full scale before the elution time of alpha-BHC, and also return to below 25% of full scale after the elution time of alpha-BHC and before the elution time of DCB?

 Yes? or No?
- 8. If a chromatogram is replotted electronically to meet these requirements; is the scaling factor used displayed on the chromatogram, and both the initial chromatogram and the replotted chromatogram submitted in the data package.

 Yes? or No?

Action:

- a. If the qualitative criteria for both columns were not met, all target compounds that are reported as detected should be considered non-detected.
- b. Use professional judgment to assign an appropriate quantitation limit using the following quidance:
 - If the detected target compound peak was sufficiently outside the pesticide RT Window, the reported values may be a false positive and should be replaced with the sample Contract Required Quantitation Limits (CRQL) value.

- ii. If the detected target compound peak poses an interference with potential detection of another target peak, the reported value should be considered and qualified as unusable (R).
- c. If the data reviewer identifies a peak in both GC column analyses that falls within the appropriate RT Windows, but was reported as a non-detect, the compound may be a false negative. Use professional judgment to decide if the compound should be included.

Note: State in the Data Review Narrative all conclusions made regarding target compound identification.

- d. If the Toxaphene peak RT windows determined from the calibration overlap with SCPs or chromatographic interferences, use professional judgment to qualify the data.
- e. If target compounds were detected on both GC columns, and the Percent Difference between the two results is greater than 25.0%, consider the potential for coelution and use professional judgment to decide whether a much larger concentration obtained on one column versus the other indicates the presence of an interfering compound. If an interfering compound is indicated, use professional judgment to determine how best to report, and if necessary, qualify the data according to these guidelines.
- f. If Toxaphene exhibits a marginal pattern-matching quality, use professional judgment to establish whether the differences are due to environmental "weathering" (i.e., degradation of the earlier eluting peaks relative to the later eluting peaks). If the presence of Toxaphene is strongly suggested, report results as presumptively present (N).

GAS CHROMATOGRAPH/MASS SPECTROMETER (GC/MS) CONFIRMATION

NOTE: This confirmation is not usually provided by the laboratory. In cases where it is provided, use professional judgment to determine if data qualified with "C" can be salvaged if it was previously qualified as unusable (R).

Action:

- a. If the quantitative criteria for both columns were met (≥ 5.0 ng/µL for SCPs and ≥ 125 ng/µL for Toxaphene), determine whether GC/MS confirmation was performed. If it was performed, qualify the data using the following guidance:
 - i. If GC/MS confirmation was not required because the quantitative criteria for both columns was not met, but it was still performed, use professional judgment when evaluating the data to decide whether the detect should be qualified with "C".
 - ii. If GC/MS confirmation was performed, but unsuccessful for a target compound detected by GC/ECD analysis, qualify those detects as "X".

All criteria were met _	_X_	_
Criteria were not met		
and/or see below		

RF = 0.744

COMPOUND QUANTITATION AND REPORTED CONTRACT REQUIRED QUANTITATION LIMITS (CRQLS)

The sample quantitation evaluation is to verify laboratory quantitation results. In the space below, please show a minimum of one sample calculation:

JC21103-1 Decachlorobiphenyl

[] = (88237732)(50)/(369.0 X 10⁶)(0.744)

= 16.07 ppb Ok

Action:

- a. If sample quantitation is different from the reported value, qualify result as unusable (R).
- b. When a sample is analyzed at more than one dilution, the lowest CRQLs are used unless a QC exceedance dictates the use of the higher CRQLs from the diluted sample.
- c. Replace concentrations that exceed the calibration range in the original analysis by crossing out the "E" and its corresponding value on the original reporting form and substituting the data from the diluted sample.
- d. Results between the MDL and CRQL should be qualified as estimated (J).
- e. Results less than the MDL should be reported at the CRQL and qualified (U). MDLs themselves are not reported.
- f. For non-aqueous samples, if the percent moisture is less than 70.0%, no qualification of the data is necessary. If the percent moisture is greater than or equal to 70.0% and less than 90.0%, qualify detects as estimated (J) and non-detects as approximated (UJ). If the percent moisture is greater than or equal to 90.0%, qualify detects as estimated (J) and non-detects as unusable (R) (see Table).

Percent Moisture Actions for Pesticide Analysis for Non-Aqueous Samples

Criteria	Action		
	Detected Associated Non-detected Associat Compounds Compounds		
% Moisture < 70.0	No qualification		
70.0 < % Moisture < 90.0	J	UJ	
% Moisture > 90.0	J R		

t sample	s which hav	re ≤ 50 % s	olids			
	<u></u>			 		_
			<u> </u>		 	

Note: If any discrepancies are found, the Region's designated representative may contact the laboratory to obtain additional information that could resolve any differences. If a discrepancy remains unresolved, the reviewer must use professional judgment to decide which value is the most accurate. Under these circumstances, the reviewer may determine that qualification of data is warranted. Note in the Data Review Narrative a description of the reasons for data qualification and the qualification that is applied to the data.

Dilution performed

SAMPLE ID	DILUTION FACTOR	REASON FOR DILUTION
·		
	•	
	-	
		-
		1

All criteria were melN/A
Criteria were not met
and/or see below

FIELD DUPLICATE PRECISION

Sample IDs:

NOTE: In the absence of QAPP guidance for validating data from field duplicates, the following action will be taken.

Matriv.

Field duplicates samples may be taken and analyzed as an indication of overall precision. These analyses measure both field and lab precision; therefore, the results may have more variability than laboratory duplicates which only laboratory performance. It is also expected that soil duplicate results will have a greater variance than water matrices due to difficulties associated with collecting identical field duplicate samples. Identify which samples within the data package are field duplicates. Estimate the relative percent difference (RPD) between the values for each compound. If large RPDs (> 50%) is observed, confirm identification of samples and note difference in the executive summary.

Campic ID	J			WIGH IX.		
COMPOUND	SQL ug/L	SAMPLE CONC.	DUPLICATE CONC.	RPD	ACTION	
No field/laboratory duplicate analyzed with this data package. LCS/LCSD % recoveries RPD used to assess precision. RPD within the required criteria of < 50 %.						

Actions:

- a. Qualify as estimated positive results (J) and nondetects (UJ) for the compound that exceeded the above criteria. For organics, only the sample and duplicate will be qualified.
- b. If an RPD cannot be calculated because one or both of the sample results is not detected, the following actions apply:
 - If one sample result is not detected and the other is greater than 5x the SQL qualify (J/UJ).
 - ii. If one sample value is not detected and the other is greater than 5x the SQL and the SQLs for the sample and duplicate are significantly different, use professional judgment to determine if qualification is appropriate.
 - iii. If one sample value is not detected and the other is less than 5x, use professional judgment to determine if qualification is appropriate.
 - iv. If both sample and duplicate results are not detected, no action is needed.

OVERALL ASSESSMENT OF DATA

Action:

- 1. Use professional judgment to determine if there is any need to qualify data which were not qualified based on the Quality Control (QC) criteria previously discussed.
- 2. Write a brief narrative to give the user an indication of the analytical limitations of the data.

The Contract Laboratory Program Project Officer (CLP PO) must be informed if any inconsistency of the data with the Sample Delivery Group (SDG) Narrative. If sufficient information on the intended use and required quality of the data is available, the reviewer should include their assessment of the usability of the data within the given context. This may be used as part of a formal Data Quality Assessment (DQA).

Overall assessment of the data: Results are valid; the data can be used for

EXECUTIVE NARRATIVE

SDG No:

JC21103

Laboratory:

Accutest, Florida

Analysis:

SW846-8015C

Number of Samples:

0

Location:

BMSMC, Building 5 Area

Humacao, PR

SUMMARY:

Eight (8) samples were analyzed for the low molecular weight alcohols (LMWAs) list following method SW846-8015C. The sample results were assessed according to USEPA data validation guidance documents in the following order of precedence: "Test Methods for Evaluating Solid Waste, Physical/Chemical Methods SW-846 (Final Update III, December 1996)," specifically for Methods 8000/8015C are utilized. The QC criteria and data validation actions listed on the data review worksheets are from the primary guidance document, unless otherwise noted.

Results are valid and can be used for decision making purposes.

Critical issues:

None

Major:

None

Minor:

None

Critical findings:

None

Major findings:

None

Minor findings:

1. MS/MSD % recoveries RPD for ethanol outside the laboratory control limits. No action taken, professional judgment. RPD was outside the laboratory control limits but within

generally acceptable control limits. No qualification made based on RPD result.

COMMENTS:

Results are valid and can be used for decision making purposes.

Reviewers Name:

Rafael infante

Chemist License 1888

Signature:

Date:

June 22, 2016

SAMPLE ORGANIC DATA SAMPLE SUMMARY

Sample ID: JC21103-1

Sample location: BMSMC Building 5 Area

Sampling date: 5/25/2016

Matrix: Groundwater

METHOD: 8015C

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Ethanol	100	ug/l	1.0	-	U	Yes
Isobutyl Alcohol	100	ug/l	1.0	-	U	Yes
Isopropyl Alcohol	100	ug/l	1.0	-	U	Yes
n-Propyl Aicohol	100	ug/l	1.0	-	U	Yes
n-Butyl Alcohol	100	ug/l	1.0	-	U	Yes
sec-Butyl Alcohol	100	ug/l	1.0	-	U	Yes
Methanol	200	ug/l	1.0	-	U	Yes

Sample ID: JC21103-2

Sample location: BMSMC Building 5 Area

Sampling date: 5/25/2016

Matrix: AQ - Equipment Blank

METHOD: 8015C

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Ethanol	100	ug/!	1.0	-	U	Yes
Isobutyl Alcohol	100	ug/l	1.0	-	U	Yes
Isopropyl Alcohol	100	ug/l	1.0	-	U	Yes
n-Propyl Alcohol	100	ug/l	1.0	-	U	Yes
n-Butyl Alcohol	100	ug/l	1.0	-	U	Yes
sec-Butyl Alcohol	100	ug/l	1.0	-	υ	Yes
Methanol	200	ug/l	1.0	-	U	Yes

Sample ID: JC21103-3

Sample location: BMSMC Building 5 Area

Sampling date: 5/26/2016

Matrix: Soil

METHOD: 8015C

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Ethanol	110	ug/kg	1.0	-	U	Yes
Isobutyl Alcohol	110	ug/kg	1.0	-	U	Yes
Isopropyl Alcohol	110	ug/kg	1.0	-	U	Yes
n-Propyl Alcohol	110	ug/kg	1.0	•	ប	Yes
n-Butyl Alcohol	110	ug/kg	1.0	-	บ	Yes
sec-Butyl Alcohol	110	ug/kg	1.0	-	บ	Yes
Methanol	220	ue/kg	1.0	-	U	Yes

Sample ID: JC21103-4

Sample location: BMSMC Building 5 Area

Sampling date: 5/26/2016

Matrix: Soil

METHOD: 8015C

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Ethanol	110	ug/kg	1.0	-	U	Yes
Isobutyl Alcohol	110	ug/kg	1.0	-	U	Yes
Isopropyl Alcohol	110	ug/kg	1.0	-	U	Yes
n-Propyl Alcohol	110	ug/kg	1.0	-	U	Yes
n-Butyl Alcohol	110	ug/kg	1.0	-	U	Yes
sec-Butył Alcohol	110	ug/kg	1.0	-	U	Yes
Methanol	220	ug/kg	1.0	-	U	Yes

Sample ID: JC21103-5

Sample location: BMSMC Building 5 Area

Sampling date: 5/26/2016

Matrix: Soil

METHOD: 8015C

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Ethanol	110	ug/kg	1.0	-	U	Yes
Isobutyl Alcohol	110	ug/kg	1.0	-	U	Yes
Isopropyl Alcohol	110	ug/kg	1.0	-	U	Yes
n-Propyl Alcohol	110	ug/kg	1.0	-	U	Yes
n-Butyl Alcohol	110	ug/kg	1.0	-	U	Yes
sec-Butyl Alcohol	110	ug/kg	1.0	-	υ	Yes
Methanol	230	ug/kg	1.0	•	U	Yes

Sample ID: JC21103-6

Sample location: BMSMC Building 5 Area

Sampling date: 5/26/2016

Matrix: Soil

METHOD: 8015C

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Ethanol	120	ug/kg	1.0		U	Yes
Isobutyl Alcohol	120	ug/kg	1.0	-	U	Yes
Isopropyl Alcohol	120	ug/kg	1.0	-	U	Yes
n-Propyi Alcohol	120	ug/kg	1.0	-	U	Yes
n-Butyl Alcohol	120	ug/kg	1.0	-	U	Yes
sec-Butyl Alcohol	120	ug/kg	1.0	-	U	Yes
Methanol	240	ug/kg	1.0	-	U	Yes

1	Project Number:	_JC21103
1	Date:	05/25-26/2016
		_05/26/2016
		2
REVIEW OF VOLATILE ORGO The following guidelines for evaluating volatile organics were or document will assist the reviewer in using professional judgme serving the needs of the data users. The sample results we guidance documents in the following order of precedence Physical/Chemical Methods SW-846 (Final Update III, Decembe utilized. The QC criteria and data validation actions listed on guidance document, unless otherwise noted. The hardcopied (laboratory name) _Accutest and the quality control and performance data summarized. The re	reated to delineate rent to make more in the assessed accorde: "Test Methods er 1996)," specifical the data review word data package modified data review	nformed decision and in better ding to USEPA data validation for Evaluating Solid Waste, ly for Methods 8000/8015C are brksheets are from the primary the received has been reviewed by for VOCs included:
Lab. Project/SDG No.:JC21103 No. of Samples:8	Sample matrix: _	SoliGroundwater
Trip blank No.: Field blank No.: Equipment blank No.: JC21103-2 Field duplicate No.: JC21103-3/JC21103-4		
X Data CompletenessX Holding TimesN/A_ GC/MS TuningN/A_ Internal Standard Performance	X Laborato X Field Dup X Calibratio	ry Control Spikes plicates ons nd Identifications nd Quantitation
Overall Comments:_Low_molecular_weight_alcoho	ls_by_SW-846_80 	015C
Definition of Qualifiers: J- Estimated results U- Compound not detected R- Rejected data UJ- Estimated nondetect Reviewer: Date:June_22,2016		

\$) ¥)			

DATA COMPLETENESS

MISSING INFORMATION	DATE LAB. CONTACTED	DATE RECEIVED

All criteria were met _	х_
Criteria were not met	
and/or see below	

HOLDING TIMES

The objective of this parameter is to ascertain the validity of the results based on the holding time of the sample from time of collection to the time of analysis.

Complete table for all samples and note the analysis and/or preservation not within criteria

SAMPLE ID	DATE SAMPLED DATE ANALYZED pH		ACTION		
All samples analy preserved.	 yzed within the red	commended method h	nolding tir	me. All samples properly	

<u>Criteria</u>

Aqueous samples – 14 days from sample collection for preserved samples (pH \leq 2, 4°C), no air bubbles. Aqueous samples – 7 days from sample collection for unpreserved samples, 4°C, no air bubbles.

Soil samples- 7 days from sample collection.

Cooler temperature (Criteria: 4 ± 2 °C): 4°C

Actions

If the VOCs vial(s) have air bubbles, estimate positive results (J) and reject nondetects (R).

If the % solids of soil samples is 10-50%, estimates positive results (J) and nondetects (UJ)

If the % solid of soil samples is < 10%, estimate positive results (J) and reject nondetects (R).

If holding times are exceeded but < 14 days beyond criteria, estimate positive results (J) and nondetects (UJ).

If holding times are exceeded but < 28 days beyond criteria, estimate positive results (J) and reject nondetects (R).

If holding times are grossly exceeded (> 28 days beyond criteria), reject all results (R).

If samples were not iced or if the ice were melted (> 10°C), estimate positive results (J) and nondetects (UJ).

All criteria were metN/A Criteria were not met see below
f the sample instrumentation is within the standard
and found to be within the specified criteria.
of sample analysis.
the associated data should be accepted, qualified

GC/MS TUNING The assessment of the tuning results is to determine if tuning QC limits __N/A_ The BFB performance results were reviewed a __N/A_ BFB tuning was performed for every 12 hours of If no, use professional judgment to determine whether or rejected. List the samples affected: If mass calibration is in error, all associated data are rejected.

All criteria were met _X
Criteria were not met
and/or see below

CALIBRATION VERIFICATION

Compliance requirements for satisfactory instrument calibration are established to ensure that the instrument is capable of producing and maintaining acceptable quantitative data.

	Date of initial calibr	ation:05/17/	16
	Dates of continuing	calibration:_05/17/16 (init	tial);_05/27/16;_05/31/16_
	Dates of final calibr	ation verification:05/27/	16;_05/31/16
	Instrument ID numb	per:GCG	H
	Matrix/Level:	Aqueous/low_	
LAB FILE ID#	CRITERIA OUT	COMPOUND	SAMPLES
	RFs, %RSD, %D, r		AFFECTED

Note: Initial, continuing, and final calibration verifications meets method specific criteria in at least one of the column, second column used for confirmation only.

Criteria

DATE

All RFs must be > 0.05 regardless of method requirements for SPCC.

All %RSD must be < 15 % regardless of method requirements for CCC.

All %Ds must be \leq 20% regardless of method requirements for CCC.

It should be noted that Region 2 SOP HW-24 does not specify criterion for the curve correlation coefficient (r). A limit for r of ≥ 0.995 has therefore been utilized as professional judgment.

Actions

If any compound has an initial RF or a continuing RF of < 0.05, estimate positive results (J) and reject nondetects (R), regardless of method requirements.

If any compound has a %RSD > 15%, estimate positive results (J) and use professional judgment to qualify nondetects.

If any compound has a %RSD > 90%, estimate positive results (J) and reject nondetects (R).

If any compound has a % D > 20%, estimate positive results (J) and reject nondetects (R).

If any compound has a % D > 20%, estimate positive results (J) and nondetects (UJ).

If any compound has a % D > 90%, estimate positive results (J) and reject nondetects (R).

If any compound has r < 0.995, estimate positive results and nondetects.

A separate worksheet should be filled for each initial curve

All criteria were met _X
Criteria were not met
and/or see below

V A. BLANK ANALYSIS RESULTS (Sections 1 & 2)

The assessment of the blank analysis results is to determine the existence and magnitude of contamination problems. The criteria for evaluation of blanks apply only to blanks associated with the samples, including trip, equipment, and laboratory blanks. If problems with any blanks exist, all data associated with the case must be carefully evaluated to determine whether or not there is an inherent variability in the data for the case, or if the problem is an isolated occurrence not affecting other data.

List the contamination in the blanks below. High and low levels blanks must be treated separately.

Laboratory blanks

DATE ANALYZED	LAB ID	LEVEL/ MATRIX	COMPOUND	CONCENTRATION UNITS
Field/ <u>Equipmen</u>				
DATE ANALYZED	LAB ID	LEVEL/ MATRIX	COMPOUND	CONCENTRATION UNITS
			nt_blankNo_field/trip_	blanks_included_in_this_data

All criteria were met _	х_
Criteria were not met	
and/or see below	

VB. BLANK ANALYSIS RESULTS (Section 3)

Blank Actions

Action Levels (ALs) should be based upon the highest concentration of contaminant determined in any blank. Do not qualify any blank with another blank. The ALs for samples which have been diluted should be corrected for the sample dilution factor and/or % moisture, where applicable. No positive sample results should be reported unless the concentration of the compound in the samples exceeds the ALs:

ALs = 10x the amount of common contaminants (methylene chloride, acetone, 2-butanone, and toluene) ALs = 5x for any other compounds

Specific actions are as follows:

If the concentration is < sample quantitation limit (SQL) and \le AL, report the compound as not detected (U) at the SQL.

If the concentration is \geq SQL but \leq AL, report the compound as not detected (U) at the reported concentration.

If the concentration is \geq SQL and > AL, report the concentration unqualified.

Notes:

High and low level blanks must be treated separately

Compounds qualified "U" for blank contamination are still considered "hits" when qualifying for calibration criteria.

CONTAMINATION SOURCE/LEVEL	COMPOUND	CONC/UNITS	AL/UNITS	SQL	AFFECTED SAMPLES
-					
		,			

All criteria were metX
Criteria were not met
and/or see below

SURROGATE SPIKE RECOVERIES

Laboratory performance of individual samples is established by evaluation of surrogate spike recoveries. All samples are spiked with surrogate compounds prior to sample analysis. The accuracy of the analysis is measured by the surrogate percent recovery. Since the effects of the sample matrix are frequently outside the control of the laboratory and may present relatively unique problems, the validation of data is frequently subjective and demands analytical experience and professional judgment. List the percent recoveries (%Rs) which do not meet the criteria for surrogate recovery. Matrix: solid/aqueous

SAMPLE ID		SURROGATE COMPOUND				ACTION	
He	xanol	DBFM	TOL-d8	BFB			
_All_surrogate_recove	ries_within_la	boratory_contro	ol_limits				
						_	
QC Limits* (Aqueous)							
, , ,		to	to	to	<u> </u>		
LL_to_UL	_69_to_121	to	to	to			
QC Limits* (Solid-Med) LL_to_UL	to	to	to	to			
1,2-DCA = 1,2-Dichloro DBFM = Dibromofluoro				Toluene-d8 mofluorober	nzene		
 If QC limits are 			nce criteria, LL = 0 – 120 % for aq		UL = upper limit. 70 – 130 % for	solid	
samples. Actions:							
QUALITY		%R < 10%	%R = 10%	5 - LL %	R > UL		
Positive result	e l	.1	1.1				

Surrogate action should be applied:

Nondetects results

If one or more surrogate in the VOC fraction is out of specification, but has a recovery of > 10%. If any one surrogate in a fraction shows < 10 % recovery.

UJ

Accept

All criteria were met _X
Criteria were not met
and/or see below

VII. A MATRIX SPIKE/MATRIX SPIKE DUPLICATE (MS/MSD)

This data is generated to determine long term precision and accuracy in the analytical method for various matrices. This data alone cannot be used to evaluate the precision and accuracy of individual samples. If any % R in the MS or MSD falls outside the designated range, the reviewer should determine if there are matrix effects, i.e. LCS data are within the QC limits but MS/MSD data are outside QC limit.

1. MS/MSD Recoveries and Precision Criteria

The laboratory should use one MS and a duplicate analysis of an unspiked field sample if target analytes are expected in the sample. If target analytes are not expected, MS/MSD should be analyzed.

List the %Rs, RPD of the compounds which do not meet the criteria.

	1103-1MS/-MSD 1103-3AMS/-MSD			_	Groundwater/ow Soil/ow
MS OR MSD	COMPOUND	% R	RPD	QC LIMITS	ACTION
	overies_and_RPD_w nt:			control_limits_ex	cept_in_the_cases_described
_MS/MSD	Ethanol		27	22	No_action

Note: No action taken, professional judgment. RPD was outside the laboratory control limits but within generally acceptable control limits. No qualification made based on RPD result.

* QC limits are laboratory in-house performance criteria, LL = lower limit, UL = upper limit.

If QC limits are not available, use limits of 70 – 130 %.

Actions:

QUALITY	%R < LL	%R > UL
Positive results	J	J
Nondetects results	R	Accept

All criteria were met _	X_
Criteria were not met	
and/or see below	

MS/MSD criteria apply only to the unspiked sample, its dilutions, and the associated MS/MSD samples:

If the % R for the affected compounds were < LL (or 70 %), qualify positive results (J) and nondetects (UJ).

If the % R for the affected compounds were > UL (or 130 %), only qualify positive results (J). If 25 % or more of all MS/MSD %R were < LL (or 70 %) or if two or more MS/MSD %Rs were < 10%, qualify all positive results (J) and reject nondetects (R).

VII. B MATRIX SPIKE/MATRIX SPIKE DUPLICATE

MS/MSD — Unspiked Compounds

It should be noted that Region 2 SOP HW-24 does not specify a MS/MSD criteria for the unspiked compounds in the sample. A %RSD of < 50% has therefore been utilized as professional judgment.

If all target analytes were spiked in the MS/MSD, this review element is not applicable.

List the %RSD of the compounds which do not meet the criteria.

Sample ID:		Matrix/Level/Unit:				
COMPOUND	SAMPLE CONC.	MS CONC.	MSD CONC.		ACTION	
<u>-</u>						

Actions:

A separate worksheet should be used for each MS/MSD pair.

^{*} If the % RSD > 50, qualify the positive result in the unspiked samples as estimated (J).

^{*} If the % RSD is not calculated (NC) due to nondetected value, use professional judgment to qualify the data.

All criteria were met _	X
Criteria were not met	
and/or see below	

VIII. LABORATORY CONTROL SAMPLE (LCS) ANALYSIS

This data is generated to determine accuracy of the analytical method for various matrices.

1. LCS Recoveries Criteria

Where LCS spiked with the same analyte at the same concentrations as the MS/MSD? Yes or No. If no make note in data review memo.

List the %R of compounds which do not meet the criteria

	LCS ID	COMPOUND	% R	QC LIMIT				
_Recoveries_within_laboratory_control_limits								
	1				_			

- * QC limits are laboratory in-house performance criteria, LL = lower limit, UL = upper limit
- * If QC limits are not available, use limits of 70 130 %.

Actions:

QUALITY	%R < LL	%R > UL
Positive results	J	J
Nondetects results	R	Accept

All analytes in the associated sample results are qualified for the following criteria.

If 25 % of the LCS recoveries were < LL (or 70 %), qualify all positive results (j) and reject nondetects (R).

If two or more LCS were below 10 %, qualify all positive results as (J) and reject nondetects (R).

2. Frequency Criteria:

Where LCS analyzed at the required frequency and for each matrix? <u>Yes</u> or No. If no, the data may be affected. Use professional judgment to determine the severity of the effect and qualify data accordingly. Discuss any actions below and list the samples affected.

		All criteria were metX Criteria were not met and/or see below
IX.	FIELD/LABORATORY DUPLICATE PRECISION	
	Sample IDs:JC21103-3/JC21103-4	Matrix:Soil

Field/laboratory duplicates samples may be taken and analyzed as an indication of overall precision. These analyses measure both field and lab precision; therefore, the results may have more variability than laboratory duplicates which only laboratory performance. It is also expected that soil duplicate results will have a greater variance than water matrices due to difficulties associated with collecting identical field duplicate samples.

The project QAPP should be reviewed for project-specific information. Suggested criteria: RPD \pm 30% for aqueous samples, RPD \pm 50 % for solid samples. If both samples and duplicate are <5 SQL, the RPD criteria is doubled.

COMPOUND	SQL	SAMPLE CONC.	DUPLICATE CONC.	RPD	ACTION	
Field/laboratory duplicate analyzed with this data package. RPD within laboratory and generally acceptable control limits.						

Actions:

Qualify as estimated positive results (J) and nondetects (UJ) for the compound that exceeded the above criteria. For organics, only the sample and duplicate will be qualified.

If an RPD cannot be calculated because one or both of the sample results is not detected, the following actions apply:

If one sample result is not detected and the other is greater than 5x the SQL qualify (J/UJ).

If one sample value is not detected and the other is greater than 5x the SQL and the SQLs for the sample and duplicate are significantly different, use professional judgment to determine if qualification is appropriate.

If one sample value is not detected and the other is less than 5x, use professional judgment to determine if qualification is appropriate.

If both sample and duplicate results are not detected, no action is needed.

All criteria were met _	_N/A
Criteria were not met	
and/or see below	

X. INTERNAL STANDARD PERFORMANCE

The assessment of the internal standard (IS) parameter is used to assist the data reviewer in determining the condition of the analytical instrumentation.

List the internal standard area of samples which do not meet the criteria.

- * Area of +100% or -50% of the IS area in the associated calibration standard.
- * Retention time (RT) within 30 seconds of the IS area in the associated calibration standard.

DATE	SAMPLE ID	IS OUT	IS AREA	ACCEPTABLE RANGE	ACTION	
2.000			***			_
	E					_
						_
						_
	K-0		Augusti	- 0.70	V-3. U.	

Actions:

1. IS actions should be applied to the compound quantitated with the out-of-control ISs

QUALITY	IS AREA < -25%	IS AREA = -25 % TO 50%	IS AREA > + 100%
Positive results	J	J	J
Nondetected results	R	UJ	ACCEPT

2. If a IS retention time varies more than 30 seconds, the chromatographic profile for that sample must be examined to determine if any false positive or negative exists. For shifts of a large magnitude, the reviewer may consider partial or total rejection of the data for the sample fraction.

All criteria were met _	Χ_	_
Criteria were not met		
and/or see below	_	

XII. SAMPLE QUANTITATION

The sample quantitation evaluation is to verify laboratory quantitation results. In the space below, please show a minimum of one sample calculation:

JC21103-1

Hexanol

RF = 67.60

[] = (261700)/(67.60)

= 3,871 ppm OK

All criteria were met __X__ Criteria were not met and/or see below ____

XII.	QUANT	TITATIC	I M	IMIT	2
All.	QUAIVI		JIW L	_IIVII I	O

A. Dilution performed

SAMPLE ID	DILUTION FACTOR	REASON FOR DILUTION
	-	
b.	1 2	
F2		
·		
		= 1

Percent Solids
List samples which have ≤ 50 % solids

Actions:

If the % solids of a soil sample is 10-50%, estimate positive results (J) and nondetects (UJ)

If the % solids of a soil sample is < 10%, estimate positive results (J) and reject nondetects (R)