

checked by N  
8/25/16

# CERTIFICATION

SDG No: MC47204 Laboratory: Accutest, Massachusetts  
Site: BMSMC, Phase 2A Release Matrix: Groundwater  
Assessment, Humacao, PR  
Humacao, PR

**SUMMARY:** Groundwater samples (Table 1) were collected on the BMSMC facility – Phase 2A Release Assessment Area. The BMSMC facility is located in Humacao, PR. Samples were taken August 4-5, 2016 and were analyzed in Accutest Laboratory of Marlborough, Massachusetts that reported the data under SDG No.: MC47204. Results were validated using the following quality control criteria of the methods employed (MAPED EPH, Massachusetts Department of Environmental Protection, 2004) and the latest validation guidelines (July, 2015) of the EPA Hazardous Waste Support Section. The analyses performed are shown in Table 1. Individual data review worksheets are enclosed for each target analyte group. The data sample organic data samples summary form shows for analytes results that were qualified.

In summary the results are valid and can be used for decision taking purposes.

Table 1. Samples analyzed and analysis performed

SAMPLE ID	SAMPLE DESCRIPTION	MATRIX	ANALYSIS PERFORMED
MC47204-1	OSMW-4D	Groundwater	Extractable TPHC Ranges
MC47204-2	OSMW-4S	Groundwater	Extractable TPHC Ranges
MC47204-2D	OSMW-4S MSD	Groundwater	Extractable TPHC Ranges
MC47204-2S	OSMW-4S MS	Groundwater	Extractable TPHC Ranges
MC47204-3	OSMW-5D	Groundwater	Extractable TPHC Ranges
MC47204-4	OSMW-5S	Groundwater	Extractable TPHC Ranges

Reviewer Name: Rafael Infante  
Chemist License 1888

Signature:

*Rafael Infante*

Date:

August 16, 2016



SGS Accutest

## Report of Analysis

Page 1 of 1


Client Sample ID:	OSMW-4D	Date Sampled:	08/04/16
Lab Sample ID:	MC47204-1	Date Received:	08/06/16
Matrix:	AQ - Ground Water	Percent Solids:	n/a
Method:	MADEP EPH REV 1.1 SW846 3510C		
Project:	BMSMC Phase 2A Release Assessment, Humacao, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	DE15161.D	1	08/09/16	TA	08/07/16	OP48360	GDE846
Run #2							

Run #	Initial Volume	Final Volume
Run #1	990 ml	2.0 ml
Run #2		

CAS No.	Compound	Result	RL	MDL	Units	Q
83-32-9	Acenaphthene	ND	5.1	1.6	ug/l	
208-96-8	Acenaphthylene	ND	5.1	0.36	ug/l	
120-12-7	Anthracene	ND	5.1	0.58	ug/l	
56-55-3	Benzo(a)anthracene	ND	5.1	0.31	ug/l	
50-32-8	Benzo(a)pyrene	ND	5.1	0.30	ug/l	
205-99-2	Benzo(b)fluoranthene	ND	5.1	0.45	ug/l	
191-24-2	Benzo(g,h,i)perylene	ND	5.1	0.37	ug/l	
207-08-9	Benzo(k)fluoranthene	ND	5.1	0.36	ug/l	
218-01-9	Chrysene	ND	5.1	0.44	ug/l	
53-70-3	Dibenz(a,h)anthracene	ND	5.1	0.39	ug/l	
206-44-0	Fluoranthene	ND	5.1	0.34	ug/l	
86-73-7	Fluorene	ND	5.1	0.40	ug/l	
193-39-5	Indeno(1,2,3-cd)pyrene	ND	5.1	0.30	ug/l	
91-57-6	2-Methylnaphthalene	ND	5.1	0.46	ug/l	
91-20-3	Naphthalene	ND	5.1	0.97	ug/l	
85-01-8	Phenanthrene	ND	5.1	0.31	ug/l	
129-00-0	Pyrene	ND	5.1	0.60	ug/l	
	C11-C22 Aromatics (Unadj.)	ND	100	29	ug/l	
	C9-C18 Aliphatics	30.5	100	17	ug/l	JB
	C19-C36 Aliphatics	46.2	100	27	ug/l	JB
	C11-C22 Aromatics	ND	100	29	ug/l	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
84-15-1	o-Terphenyl	89%		40-140%
321-60-8	2-Fluorobiphenyl	88%		40-140%
3386-33-2	1-Chlorooctadecane	74%		40-140%
580-13-2	2-Bromonaphthalene	100%		40-140%

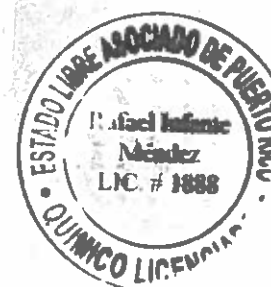

ND = Not detected      MDL = Method Detection Limit  
 RL = Reporting Limit  
 E = Indicates value exceeds calibration range

J = Indicates an estimated value  
 B = Indicates analyte found in associated method blank  
 N = Indicates presumptive evidence of a compound

SGS Accutest

## Report of Analysis

Page 1 of 1

Client Sample ID: OSMW-4S  
 Lab Sample ID: MC47204-2  
 Matrix: AQ - Ground Water  
 Method: MADEP EPH REV 1.1 SW846 3510C  
 Project: BMSMC Phase 2A Release Assessment, Humacao, PR

Date Sampled: 08/04/16  
 Date Received: 08/06/16  
 Percent Solids: n/a

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	DE15162.D	1	08/09/16	TA	08/07/16	OP48360	GDE846
Run #2							

Run #	Initial Volume	Final Volume
Run #1	975 ml	2.0 ml
Run #2		

CAS No.	Compound	Result	RL	MDL	Units	Q
83-32-9	Acenaphthene	ND	5.1	1.6	ug/l	
208-96-8	Acenaphthylene	ND	5.1	0.36	ug/l	
120-12-7	Anthracene	ND	5.1	0.59	ug/l	
56-55-3	Benzo(a)anthracene	ND	5.1	0.31	ug/l	
50-32-8	Benzo(a)pyrene	ND	5.1	0.30	ug/l	
205-99-2	Benzo(b)fluoranthene	ND	5.1	0.46	ug/l	
191-24-2	Benzo(g,h,i)perylene	ND	5.1	0.38	ug/l	
207-08-9	Benzo(k)fluoranthene	ND	5.1	0.36	ug/l	
218-01-9	Chrysene	ND	5.1	0.44	ug/l	
53-70-3	Dibenz(a,h)anthracene	ND	5.1	0.40	ug/l	
206-44-0	Fluoranthene	ND	5.1	0.34	ug/l	
86-73-7	Fluorene	ND	5.1	0.41	ug/l	
193-39-5	Indeno(1,2,3-cd)pyrene	ND	5.1	0.30	ug/l	
91-57-6	2-Methylnaphthalene	ND	5.1	0.46	ug/l	
91-20-3	Naphthalene	ND	5.1	0.98	ug/l	
85-01-8	Phenanthrene	ND	5.1	0.31	ug/l	
129-00-0	Pyrene	ND	5.1	0.61	ug/l	
	C11-C22 Aromatics (Unadj.)	ND	100	29	ug/l	
	C9-C18 Aliphatics	34.4	100	17	ug/l	JB
	C19-C36 Aliphatics	49.4	100	28	ug/l	JB
	C11-C22 Aromatics	ND	100	29	ug/l	

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
84-15-1	o-Terphenyl	70%		40-140%
321-60-8	2-Fluorobiphenyl	81%		40-140%
3386-33-2	1-Chlorooctadecane	68%		40-140%
580-13-2	2-Bromonaphthalene	92%		40-140%



ND = Not detected      MDL = Method Detection Limit  
 RL = Reporting Limit  
 E = Indicates value exceeds calibration range

J = Indicates an estimated value  
 B = Indicates analyte found in associated method blank  
 N = Indicates presumptive evidence of a compound

SGS Accutest

## Report of Analysis

Page 1 of 1

Client Sample ID:	OSMW-5D	Date Sampled:	08/05/16
Lab Sample ID:	MC47204-3	Date Received:	08/06/16
Matrix:	AQ - Ground Water	Percent Solids:	n/a
Method:	MADEP EPH REV 1.1 SW846 3510C		
Project:	BMSMC Phase 2A Release Assessment, Humacao, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	DE15163.D	1	08/09/16	TA	08/07/16	OP48360	GDE846
Run #2							

Run #	Initial Volume	Final Volume
Run #1	960 ml	2.0 ml
Run #2		

CAS No.	Compound	Result	RL	MDL	Units	Q
83-32-9	Acenaphthene	ND	5.2	1.6	ug/l	
208-96-8	Acenaphthylene	ND	5.2	0.37	ug/l	
120-12-7	Anthracene	ND	5.2	0.60	ug/l	
56-55-3	Benzo(a)anthracene	ND	5.2	0.32	ug/l	
50-32-8	Benzo(a)pyrene	ND	5.2	0.30	ug/l	
205-99-2	Benzo(b)fluoranthene	ND	5.2	0.47	ug/l	
191-24-2	Benzo(g,h,i)perylene	ND	5.2	0.39	ug/l	
207-08-9	Benzo(k)fluoranthene	ND	5.2	0.37	ug/l	
218-01-9	Chrysene	ND	5.2	0.45	ug/l	
53-70-3	Dibenz(a,h)anthracene	ND	5.2	0.40	ug/l	
206-44-0	Fluoranthene	ND	5.2	0.35	ug/l	
86-73-7	Fluorene	ND	5.2	0.41	ug/l	
193-39-5	Indeno(1,2,3-cd)pyrene	ND	5.2	0.30	ug/l	
91-57-6	2-Methylnaphthalene	ND	5.2	0.47	ug/l	
91-20-3	Naphthalene	ND	5.2	1.0	ug/l	
85-01-8	Phenanthrene	ND	5.2	0.32	ug/l	
129-00-0	Pyrene	ND	5.2	0.62	ug/l	
	C11-C22 Aromatics (Unadj.)	36.0	100	30	ug/l	J
	C9-C18 Aliphatics	38.8	100	17	ug/l	JB
	C19-C36 Aliphatics	51.2	100	28	ug/l	JB
	C11-C22 Aromatics	36.0	100	30	ug/l	J

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
84-15-1	o-Terphenyl	82%		40-140%
321-60-8	2-Fluorobiphenyl	81%		40-140%
3386-33-2	1-Chlorooctadecane	89%		40-140%
580-13-2	2-Bromonaphthalene	76%		40-140%



ND = Not detected      MDL = Method Detection Limit  
 RL = Reporting Limit  
 E = Indicates value exceeds calibration range

J = Indicates an estimated value  
 B = Indicates analyte found in associated method blank  
 N = Indicates presumptive evidence of a compound

SGS Accutest

## Report of Analysis

Page 1 of 1

Client Sample ID:	OSMW-5S	Date Sampled:	08/05/16
Lab Sample ID:	MC47204-4	Date Received:	08/06/16
Matrix:	AQ - Ground Water	Percent Solids:	n/a
Method:	MADEP EPH REV 1.1 SW846 3510C		
Project:	BMSMC Phase 2A Release Assessment, Humacao, PR		

Run #	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
Run #1	DE15164.D	1	08/09/16	TA	08/07/16	OP48360	GDE846
Run #2							

Run #	Initial Volume	Final Volume
Run #1	980 ml	2.0 ml
Run #2		

CAS No.	Compound	Result	RL	MDL	Units	Q
83-32-9	Acenaphthene	ND	5.1	1.6	ug/l	
208-96-8	Acenaphthylene	ND	5.1	0.36	ug/l	
120-12-7	Anthracene	ND	5.1	0.59	ug/l	
56-55-3	Benzo(a)anthracene	ND	5.1	0.31	ug/l	
50-32-8	Benzo(a)pyrene	ND	5.1	0.30	ug/l	
205-99-2	Benzo(b)fluoranthene	ND	5.1	0.46	ug/l	
191-24-2	Benzo(g,h,i)perylene	ND	5.1	0.38	ug/l	
207-08-9	Benzo(k)fluoranthene	ND	5.1	0.36	ug/l	
218-01-9	Chrysene	ND	5.1	0.44	ug/l	
53-70-3	Dibenz(a,h)anthracene	ND	5.1	0.40	ug/l	
206-44-0	Fluoranthene	ND	5.1	0.34	ug/l	
86-73-7	Fluorene	ND	5.1	0.40	ug/l	
193-39-5	Indeno(1,2,3-cd)pyrene	ND	5.1	0.30	ug/l	
91-57-6	2-Methylnaphthalene	ND	5.1	0.46	ug/l	
91-20-3	Naphthalene	ND	5.1	0.98	ug/l	
85-01-8	Phenanthrene	ND	5.1	0.31	ug/l	
129-00-0	Pyrene	ND	5.1	0.61	ug/l	
	C11-C22 Aromatics (Unadj.)	36.2	100	29	ug/l	J
	C9-C18 Aliphatics	41.7	100	17	ug/l	JB
	C19-C36 Aliphatics	55.1	100	28	ug/l	JB
	C11-C22 Aromatics	36.2	100	29	ug/l	J

CAS No.	Surrogate Recoveries	Run# 1	Run# 2	Limits
84-15-1	o-Terphenyl	86%		40-140%
321-60-8	2-Fluorobiphenyl	83%		40-140%
3386-33-2	1-Chlorooctadecane	78%		40-140%
580-13-2	2-Bromonaphthalene	80%		40-140%



ND = Not detected      MDL = Method Detection Limit  
 RL = Reporting Limit  
 E = Indicates value exceeds calibration range

J = Indicates an estimated value  
 B = Indicates analyte found in associated method blank  
 N = Indicates presumptive evidence of a compound

## Matrix Spike/Matrix Spike Duplicate Summary

Page 1 of 1

Job Number: MC47204

Account: AMANYWP Anderson Mulholland and Assoc.

Project: BMSMC Phase 2A Release Assessment, Humacao, PR

Sample	File ID	DF	Analyzed	By	Prep Date	Prep Batch	Analytical Batch
OP48360-MS	DE15159.D	1	08/09/16	TA	08/07/16	OP48360	GDE846
OP48360-MSD	DE15160.D	1	08/09/16	TA	08/07/16	OP48360	GDE846
MC47204-2	DE15162.D	1	08/09/16	TA	08/07/16	OP48360	GDE846

The QC reported here applies to the following samples:

Method: MADEP EPH REV 1.1

MC47204-1, MC47204-2, MC47204-3, MC47204-4

CAS No.	Compound	MC47204-2 ug/l	Spike Q	Spike ug/l	MS ug/l	MS %	Spike ug/l	MSD ug/l	MSD %	RPD	Limits Rec/RPD
83-32-9	Acenaphthene	ND		52.6	40.7	77	51.5	35.9	70	13	40-140/25
208-96-8	Acenaphthylene	ND		52.6	37.0	70	51.5	33.7	65	9	40-140/25
120-12-7	Anthracene	ND		52.6	38.4	73	51.5	34.3	67	11	40-140/25
56-55-3	Benzo(a)anthracene	ND		52.6	47.7	91	51.5	43.0	83	10	40-140/25
50-32-8	Benzo(a)pyrene	ND		52.6	46.2	88	51.5	41.9	81	10	40-140/25
205-99-2	Benzo(b)fluoranthene	ND		52.6	47.3	90	51.5	43.1	84	9	40-140/25
191-24-2	Benzo(g,h,i)perylene	ND		52.6	51.3	97	51.5	45.4	88	12	40-140/25
207-08-9	Benzo(k)fluoranthene	ND		52.6	45.9	87	51.5	41.5	81	10	40-140/25
218-01-9	Chrysene	ND		52.6	44.7	85	51.5	40.6	79	10	40-140/25
53-70-3	Dibenz(a,h)anthracene	ND		52.6	49.6	94	51.5	44.5	86	11	40-140/25
206-44-0	Fluoranthene	ND		52.6	45.8	87	51.5	41.5	81	10	40-140/25
86-73-7	Fluorene	ND		52.6	40.7	77	51.5	36.6	71	11	40-140/25
193-39-5	Indeno(1,2,3-cd)pyrene	ND		52.6	47.6	90	51.5	43.1	84	10	40-140/25
91-57-6	2-Methylnaphthalene	ND		52.6	39.1	74	51.5	35.8	69	9	40-140/25
91-20-3	Naphthalene	ND		52.6	32.0	61	51.5	29.4	57	8	40-140/25
85-01-8	Phenanthrene	ND		52.6	43.8	83	51.5	39.1	76	11	40-140/25
129-00-0	Pyrene	ND		52.6	44.9	85	51.5	40.8	79	10	40-140/25
	C11-C22 Aromatics (Unadj.)	ND		842	800	95	825	726	88	10	40-140/25
	C9-C18 Aliphatics	34.4	JB	316	301	84	309	258	72	15	40-140/25
	C19-C36 Aliphatics	49.4	JB	421	447	94	412	396	84	12	40-140/25

CAS No.	Surrogate Recoveries	MS	MSD	MC47204-2	Limits
84-15-1	o-Terphenyl	92%	83%	70%	40-140%
321-60-8	2-Fluorobiphenyl	88%	89%	81%	40-140%
3386-33-2	1-Chlorooctadecane	91%	80%	68%	40-140%
580-13-2	2-Bromonaphthalene	92%	89%	92%	40-140%



\* = Outside of Control Limits.



## EXECUTIVE NARRATIVE

SDG No: **MC47204** Laboratory: **Accutest, Massachusetts**  
Analysis: **MADEP EPH** Number of Samples: **6**  
Location: **BMSMC, Phase 2A Release Assessment Area**  
**Humacao, PR**

**SUMMARY:** Six (6) samples were analyzed for Volatiles TPHC Ranges by method MADEP EPH. Samples were validated following the METHOD FOR THE DETERMINATION OF EXTRACTABLE PETROLEUM HYDROCARBONS (EPH) quality control criteria, Massachusetts Department of Environmental Protection, Revision 1.1 (2004). Also the general validation guidelines promulgated by the USEPA Hazardous Wastes Support Section. The QC criteria and data validation actions listed on the data review worksheets are from the primary guidance document, unless otherwise noted.

Results are valid and can be used for decision making purposes.

**Critical issues:** **None**  
**Major:** **None**  
**Minor:** **None**

**Critical findings:** **None**  
**Major findings:** **None**  
**Minor findings:** 1. C9 – C18 aliphatics and C19 – C36 aliphatics found in method blank. No action taken, analyte concentration below reporting limit. Target analytes not detected above reporting limit in sample batch. Laboratory qualified the results with a B. No further qualification performed

**COMMENTS:** Results are valid and can be used for decision making purposes.

**Reviewers Name:** **Rafael Infante**  
**Chemist License 1888**

**Signature:**

A handwritten signature in blue ink, reading "Rafael Infante", is written over a horizontal line.

**Date:**

**August 16, 2016**



# **SAMPLE ORGANIC DATA SAMPLE SUMMARY**

**Sample ID: MC47204-1**

**Sample location: BMSMC Phase 2A Release Assessment, Humacao, PR**

**Sampling date: 8/4/2016**

**Matrix: Groundwater**

## **METHOD: 8270D**

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Acenaphthene	5.1	ug/l	1	-	U	Yes
Acenaphthylene	5.1	ug/l	1	-	U	Yes
Anthracene	5.1	ug/l	1	-	U	Yes
Benzo(a)anthracene	5.1	ug/l	1	-	U	Yes
Benzo(a)pyrene	5.1	ug/l	1	-	U	Yes
Benzo(b)fluoranthene	5.1	ug/l	1	-	U	Yes
Benzo(g,h,i)perylene	5.1	ug/l	1	-	U	Yes
Benzo(k)fluoranthene	5.1	ug/l	1	-	U	Yes
Chrysene	5.1	ug/l	1	-	U	Yes
Dibenzo(a,h)anthracene	5.1	ug/l	1	-	U	Yes
Fluoranthene	5.1	ug/l	1	-	U	Yes
Fluorene	5.1	ug/l	1	-	U	Yes
Indeno(1,2,3-cd)pyrene	5.1	ug/l	1	-	U	Yes
2-Methylnaphthalene	5.1	ug/l	1	-	U	Yes
Naphthalene	5.1	ug/l	1	-	U	Yes
Phenanthrene	5.1	ug/l	1	-	U	Yes
Pyrene	5.1	ug/l	1	-	U	Yes
C11-C22 Aromatics (Unadj.)	100	ug/l	1	-	U	Yes
C9-C18 Aliphatics	30.5	ug/l	1	JB	JB	Yes
C19-C36 Aliphatics	46.2	ug/l	1	JB	JB	Yes
C11-C22 Aromatics	100	ug/l	1	-	U	Yes

Sample ID: MC47204-2  
 Sample location: BMSMC Phase 2A Release Assessment, Humacao, PR  
 Sampling date: 8/4/2016  
 Matrix: Groundwater

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Acenaphthene	5.1	ug/l	1	-	U	Yes
Acenaphthylene	5.1	ug/l	1	-	U	Yes
Anthracene	5.1	ug/l	1	-	U	Yes
Benzo(a)anthracene	5.1	ug/l	1	-	U	Yes
Benzo(a)pyrene	5.1	ug/l	1	-	U	Yes
Benzo(b)fluoranthene	5.1	ug/l	1	-	U	Yes
Benzo(g,h,i)perylene	5.1	ug/l	1	-	U	Yes
Benzo(k)fluoranthene	5.1	ug/l	1	-	U	Yes
Chrysene	5.1	ug/l	1	-	U	Yes
Dibenzo(a,h)anthracene	5.1	ug/l	1	-	U	Yes
Fluoranthene	5.1	ug/l	1	-	U	Yes
Fluorene	5.1	ug/l	1	-	U	Yes
Indeno(1,2,3-cd)pyrene	5.1	ug/l	1	-	U	Yes
2-Methylnaphthalene	5.1	ug/l	1	-	U	Yes
Naphthalene	5.1	ug/l	1	-	U	Yes
Phenanthrene	5.1	ug/l	1	-	U	Yes
Pyrene	5.1	ug/l	1	-	U	Yes
C11-C22 Aromatics (Unadj.)	100	ug/l	1	-	U	Yes
C9-C18 Aliphatics	34.4	ug/l	1	JB	JB	Yes
C19-C36 Aliphatics	49.4	ug/l	1	JB	JB	Yes
C11-C22 Aromatics	100	ug/l	1	-	U	Yes

Sample ID: MC47204-3  
Sample location: BMSMC Phase 2A Release Assessment, Humacao, PR  
Sampling date: 8/5/2016  
Matrix: Groundwater

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Acenaphthene	5.2	ug/l	1	-	U	Yes
Acenaphthylene	5.2	ug/l	1	-	U	Yes
Anthracene	5.2	ug/l	1	-	U	Yes
Benzo(a)anthracene	5.2	ug/l	1	-	U	Yes
Benzo(a)pyrene	5.2	ug/l	1	-	U	Yes
Benzo(b)fluoranthene	5.2	ug/l	1	-	U	Yes
Benzo(g,h,i)perylene	5.2	ug/l	1	-	U	Yes
Benzo(k)fluoranthene	5.2	ug/l	1	-	U	Yes
Chrysene	5.2	ug/l	1	-	U	Yes
Dibenzo(a,h)anthracene	5.2	ug/l	1	-	U	Yes
Fluoranthene	5.2	ug/l	1	-	U	Yes
Fluorene	5.2	ug/l	1	-	U	Yes
Indeno(1,2,3-cd)pyrene	5.2	ug/l	1	-	U	Yes
2-Methylnaphthalene	5.2	ug/l	1	-	U	Yes
Naphthalene	5.2	ug/l	1	-	U	Yes
Phenanthrene	5.2	ug/l	1	-	U	Yes
Pyrene	5.2	ug/l	1	-	U	Yes
C11-C22 Aromatics (Unadj.)	36.0	ug/l	1	J	J	Yes
C9-C18 Aliphatics	38.8	ug/l	1	JB	JB	Yes
C19-C36 Aliphatics	51.2	ug/l	1	JB	JB	Yes
C11-C22 Aromatics	36.0	ug/l	1	J	J	Yes

Sample ID: MC47204-4  
 Sample location: BMSMC Phase 2A Release Assessment, Humacao, PR  
 Sampling date: 8/5/2016  
 Matrix: Groundwater

**METHOD: 8270D**

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Acenaphthene	5.1	ug/l	1	-	U	Yes
Acenaphthylene	5.1	ug/l	1	-	U	Yes
Anthracene	5.1	ug/l	1	-	U	Yes
Benzo(a)anthracene	5.1	ug/l	1	-	U	Yes
Benzo(a)pyrene	5.1	ug/l	1	-	U	Yes
Benzo(b)fluoranthene	5.1	ug/l	1	-	U	Yes
Benzo(g,h,i)perylene	5.1	ug/l	1	-	U	Yes
Benzo(k)fluoranthene	5.1	ug/l	1	-	U	Yes
Chrysene	5.1	ug/l	1	-	U	Yes
Dibenzo(a,h)anthracene	5.1	ug/l	1	-	U	Yes
Fluoranthene	5.1	ug/l	1	-	U	Yes
Fluorene	5.1	ug/l	1	-	U	Yes
Indeno(1,2,3-cd)pyrene	5.1	ug/l	1	-	U	Yes
2-Methylnaphthalene	5.1	ug/l	1	-	U	Yes
Naphthalene	5.1	ug/l	1	-	U	Yes
Phenanthrene	5.1	ug/l	1	-	U	Yes
Pyrene	5.1	ug/l	1	-	U	Yes
C11-C22 Aromatics (Unadj.)	36.2	ug/l	1	J	J	Yes
C9-C18 Aliphatics	41.7	ug/l	1	JB	JB	Yes
C19-C36 Aliphatics	55.1	ug/l	1	JB	JB	Yes
C11-C22 Aromatics	36.2	ug/l	1	J	J	Yes

Sample ID: MC47204-2MS  
Sample location: BSMC Phase 2A Release Assessment, Humacao, PR  
Sampling date: 8/4/2016  
Matrix: Groundwater

**METHOD: 8270D**

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Acenaphthene	40.7	ug/l	1	-	-	Yes
Acenaphthylene	37.0	ug/l	1	-	-	Yes
Anthracene	38.4	ug/l	1	-	-	Yes
Benzo(a)anthracene	47.7	ug/l	1	-	-	Yes
Benzo(a)pyrene	46.2	ug/l	1	-	-	Yes
Benzo(b)fluoranthene	47.3	ug/l	1	-	-	Yes
Benzo(g,h,i)perylene	51.3	ug/l	1	-	-	Yes
Benzo(k)fluoranthene	45.9	ug/l	1	-	-	Yes
Chrysene	44.7	ug/l	1	-	-	Yes
Dibenzo(a,h)anthracene	49.6	ug/l	1	-	-	Yes
Fluoranthene	45.8	ug/l	1	-	-	Yes
Fluorene	40.7	ug/l	1	-	-	Yes
Indeno(1,2,3-cd)pyrene	47.6	ug/l	1	-	-	Yes
2-Methylnaphthalene	39.1	ug/l	1	-	-	Yes
Naphthalene	32.0	ug/l	1	-	-	Yes
Phenanthrene	43.8	ug/l	1	-	-	Yes
Pyrene	44.9	ug/l	1	-	-	Yes
C11-C22 Aromatics (Unadj.)	800	ug/l	1	-	-	Yes
C9-C18 Aliphatics	301	ug/l	1	-	-	Yes
C19-C36 Aliphatics	421	ug/l	1	-	-	Yes

Sample ID: MC47204-2MSD  
Sample location: BMSMC Phase 2A Release Assessment, Humacao, PR  
Sampling date: 8/4/2016  
Matrix: Groundwater

METHOD: 8270D

Analyte Name	Result	Units	Dilution Factor	Lab Flag	Validation	Reportable
Acenaphthene	35.9	ug/l	1	-	-	Yes
Acenaphthylene	33.7	ug/l	1	-	-	Yes
Anthracene	34.3	ug/l	1	-	-	Yes
Benzo(a)anthracene	43.0	ug/l	1	-	-	Yes
Benzo(a)pyrene	41.9	ug/l	1	-	-	Yes
Benzo(b)fluoranthene	43.1	ug/l	1	-	-	Yes
Benzo(g,h,i)perylene	45.4	ug/l	1	-	-	Yes
Benzo(k)fluoranthene	41.5	ug/l	1	-	-	Yes
Chrysene	40.6	ug/l	1	-	-	Yes
Dibenzo(a,h)anthracene	44.5	ug/l	1	-	-	Yes
Fluoranthene	41.5	ug/l	1	-	-	Yes
Fluorene	36.6	ug/l	1	-	-	Yes
Indeno(1,2,3-cd)pyrene	43.1	ug/l	1	-	-	Yes
2-Methylnaphthalene	35.8	ug/l	1	-	-	Yes
Naphthalene	29.4	ug/l	1	-	-	Yes
Phenanthrene	39.1	ug/l	1	-	-	Yes
Pyrene	40.8	ug/l	1	-	-	Yes
C11-C22 Aromatics (Unadj.)	726	ug/l	1	-	-	Yes
C9-C18 Aliphatics	258	ug/l	1	-	-	Yes
C19-C36 Aliphatics	396	ug/l	1	-	-	Yes

# DATA REVIEW WORKSHEETS

Type of validation Full: ☒ Limited: ☐ Project Number: MC47204 Date: 8/04-05/2016 Shipping date: 08/05/2016 EPA Region: 2

## REVIEW OF EXTRACTABLE PETROLEUM HYDROCARBON (EPHs) PACKAGE

The following guidelines for evaluating volatile organics were created to delineate required validation actions. This document will assist the reviewer in using professional judgment to make more informed decision and in better serving the needs of the data users. The sample results were assessed according to the data validation guidance documents in the following order of precedence METHOD FOR THE DETERMINATION OF EXTRACTABLE PETROLEUM HYDROCARBONS (VPH), Massachusetts Department of Environmental Protection, Revision 1.1 (2004). Also the general validation guidelines promulgated by the USEPA Hazardous Wastes Support Section. The QC criteria and data validation actions listed on the data review worksheets are from the primary guidance document, unless otherwise noted.

The hardcopied (laboratory name) Accutest Laboratories data package received has been reviewed and the quality control and performance data summarized. The data review for SVOCs included:

Lab. Project/SDG No.: MC47204 Sample matrix: Groundwater  
No. of Samples: 6  
Field blank No.:  
Equipment blank No.:  
Trip blank No.:  
Field duplicate No.:

<input checked="" type="checkbox"/> Data Completeness	<input checked="" type="checkbox"/> Laboratory Control Spikes
<input checked="" type="checkbox"/> Holding Times	<input checked="" type="checkbox"/> Field Duplicates
<input type="checkbox"/> GC/MS Tuning	<input checked="" type="checkbox"/> Calibrations
<input type="checkbox"/> Internal Standard Performance	<input checked="" type="checkbox"/> Compound Identifications
<input checked="" type="checkbox"/> Blanks	<input checked="" type="checkbox"/> Compound Quantitation
<input checked="" type="checkbox"/> Surrogate Recoveries	<input checked="" type="checkbox"/> Quantitation Limits
<input checked="" type="checkbox"/> Matrix Spike/Matrix Spike Duplicate	

Overall Comments:  
Extractable Petroleum Hydrocarbons by GC by Method MADEP\_EPH\_REV\_1.1  
(C9\_to\_C36\_Aliphatics; C11\_to\_C22\_(Aromatics))

### Definition of Qualifiers:

J- Estimated results  
U- Compound not detected  
R- Rejected data  
UJ- Estimated nondetect

Reviewer:   
Date: 08/16/2016

## DATA REVIEW WORKSHEETS

All criteria were met   x    
Criteria were not met and/or see below           

## I. DATA COMPLETNESS

**A. Data Package:**

### MISSING INFORMATION

DATE LAB. CONTACTED

DATE RECEIVED

[illegible]

**B. Other**

**Discrepancies:**

[illegible]



## DATA REVIEW WORKSHEETS

All criteria were met   X    
Criteria were not met and/or see below           

### HOLDING TIMES

The objective of this parameter is to ascertain the validity of the results based on the holding time of the sample from time of collection to the time of extraction, and subsequently from the time of extraction to the time of analysis.

Complete table for all samples and note the analysis and/or preservation not within criteria

SAMPLE ID	DATE SAMPLED	DATE EXTRACTED	DATE ANALYZED	ACTION
Samples extracted and analyzed within method recommended holding time				

### Criteria

#### Preservation:

Aqueous samples must be acidified to a pH of 2.0 or less at the time of collection.

Soil samples must be cooled at  $4 \pm 2$  °C immediately after collection.

#### Holding times:

Samples must be extracted within 14 days of collection, and analyzed within 40 days of extraction.

Cooler temperature (Criteria:  $4 \pm 2$  °C):   2.1°C  

#### Actions: Qualify positive results/nondetects as follows:

If holding times are exceeded, estimate positive results (J) and nondetects (UJ).

If holding times are grossly exceeded, use professional judgment to qualify data. The data reviewer may choose to estimate positive results (J) and rejects nondetects (R).

If samples were not at the proper temperature ( $> 10^{\circ}\text{C}$ ) or improperly preserved, use professional judgment to qualify the results.

## DATA REVIEW WORKSHEETS

All criteria were met   X    
Criteria were not met and/or see below           

### CALIBRATIONS VERIFICATION

Compliance requirements for satisfactory instrument calibration are established to ensure that the instrument is capable of producing and maintaining acceptable quantitative data.

Date of initial calibration:           08/05/16          

Dates of initial calibration verification:           08/05/13          

Instrument ID numbers:           GCDE          

Matrix/Level:           AQUEOUS/MEDIUM          

DATE	LAB FILE ID#	ANALYTE	CRITERIA OUT RFs, %RSD, %D, r	SAMPLES AFFECTED
Initial and continuing calibration meet method specific requirements				

#### Criteria- ICAL

- Five point calibration curve.
- The percent relative standard deviation (%RSD) of the calibration factor must be equal to or less than 25% over the working range for the analyte of interest. When this condition is met, linearity through the origin may be assumed, and the average calibration factor is used in lieu of a calibration curve.
- A collective calibration factor must also be established for each hydrocarbon range of interest. Calculate the collective CFs for C9-C18 Aliphatic Hydrocarbons, C19-C36 Aliphatic Hydrocarbons, and C11-C22 Aromatic Hydrocarbons using the FID chromatogram. Tabulate the summation of the peak areas of all components in that fraction against the total concentration injected. The %RSD of the calibration factor must be equal to or less than 25% over the working range for the hydrocarbon range of interest.
  - The area for the surrogates must be subtracted from the area summation of the range in which they elute.
  - The areas associated with naphthalene and 2-methylnaphthalene in the aliphatic range standard must be subtracted from the uncorrected collective C9-C18 Aliphatic Hydrocarbon range area prior to calculating the CF.

## DATA REVIEW WORKSHEETS

### Criteria- CCAL

- At a minimum, the working calibration factor must be verified on each working day, after every 20 samples or every 24 hours (whichever is more frequent), and at the end of the analytical sequence by the injection of a mid-level continuing calibration standard to verify instrument performance and linearity.
- If the percent difference (%D) for any analyte varies from the predicted response by more than  $\pm 25\%$ , a new five-point calibration must be performed for that analyte. Greater percent differences are permissible for n-nonane. If the %D for n-nonane is greater than 30, note the nonconformance in the case narrative. It should be noted that the %Ds are calculated when CFs are used for the initial calibration and percent drifts are calculated when calibration curves using linear regression are used for the initial calibration.

### Actions:

If %RSD > 25% for target compounds or a correlation coefficient < 0.99, estimate positive results (J) and use professional judgment to qualify nondetects.

If % D > 25% (> 30 for nonane), estimate positive results (J) and nondetects (UJ).

### CALIBRATIONS VERIFICATION

Compliance requirements for satisfactory instrument calibration are established to ensure that the instrument is capable of producing and maintaining acceptable quantitative data.

Date of initial calibration: \_\_\_\_\_ 08/05/16 \_\_\_\_\_

Dates of continuing calibration verification: \_\_\_\_\_ 08/09/15 \_\_\_\_\_

Dates of final calibration verification: \_\_\_\_\_ 08/09/16 \_\_\_\_\_

Instrument ID numbers: \_\_\_\_\_ GCDE \_\_\_\_\_

Matrix/Level: \_\_\_\_\_ AQUEOUS/MEDIUM \_\_\_\_\_

DATE	LAB FILE ID#	ANALYTE	CRITERIA OUT RFs, %RSD, %D, r	SAMPLES AFFECTED
Initial and continuing calibration meet method specific requirements				

A separate worksheet should be filled for each initial curve

## DATA REVIEW WORKSHEETS

All criteria were met   X    
Criteria were not met and/or see below           

#### V A. BLANK ANALYSIS RESULTS (Sections 1 & 2)

The assessment of the blank analysis results is to determine the existence and magnitude of contamination problems. The criteria for evaluation of blanks apply only to blanks associated with the samples, including trip, equipment, and laboratory blanks. If problems with any blanks exist, all data associated with the case must be carefully evaluated to determine whether or not there is an inherent variability in the data for the case, or if the problem is an isolated occurrence not affecting other data. A Laboratory Method Blank must be run after samples suspected of being highly contaminated to determine if sample carryover has occurred.

List the contamination in the blanks below. High and low levels blanks must be treated separately.

### Laboratory blanks

DATE ANALYZED	LAB ID	LEVEL/ MATRIX	COMPOUND	CONCENTRATION UNITS
---------------	--------	---------------	----------	---------------------

METHOD BLANKS MEET THE METHOD SPECIFIC CRITERIA EXCEPT IN THE CASES DESCRIBED IN THIS DOCUMENT.

08/09/16	OP48360-MB	Aq./low	C9-C18 Aliphatics	34.8 ug/l
			C19-C36 Aliphatics	55.9 ug/l

**Note:** No action taken, analyte concentration below reporting limits. Target analyte not detected above reporting limit in sample batch. Laboratory qualified the results with a B. No further qualification performed.

### Field/Trip/Equipment

DATE ANALYZED	LAB ID	LEVEL/ MATRIX	COMPOUND	CONCENTRATION UNITS
---------------	--------	---------------	----------	---------------------

NO TRIP/FIELD/EQUIPMENT BLANKS ANALYZED ASSOCIATED WITH THIS DATA PACKAGE.

**Note:**

## DATA REVIEW WORKSHEETS

All criteria were met   X    
Criteria were not met and/or see below           

### V      B.      BLANK ANALYSIS RESULTS (Section 3)

#### Blank Actions

The ALs for samples which have been diluted should be corrected for the sample dilution factor and/or % moisture, where applicable. Peaks must not be detected above the Reporting Limit within the retention time window of any analyte of interest. The hydrocarbon ranges must not be detected at a concentration greater than 10% of the most stringent MCP cleanup standard. Specific actions area as follows:

If the concentration is < sample quantitation limit (SQL) and < AL, report the compound as not detected (U) at the SQL.

If the concentration is  $\geq$  SQL but < AL, report the compound as not detected (U) at the reported concentration.

If the concentration is > AL, report the concentration unqualified.

## DATA REVIEW WORKSHEETS

All criteria were met   X    
Criteria were not met and/or see below

## SURROGATE SPIKE RECOVERIES

Laboratory performance of individual samples is established by evaluation of surrogate spike recoveries. All samples are spiked with surrogate compounds prior to sample analysis. The accuracy of the analysis is measured by the surrogate percent recovery. Since the effects of the sample matrix are frequently outside the control of the laboratory and may present relatively unique problems, the validation of data is frequently subjective and demands analytical experience and professional judgment.

List the percent recoveries (%Rs) which do not meet the criteria for surrogate recovery.

**Matrix:** solid/aqueous

SAMPLE ID	SURROGATE COMPOUND				ACTION
	S1	S2	S3	S4	

\_SURROGATE\_STANDARDS\_RECOVERIES\_WITHIN\_LABORATORY\_CONTROL\_  
\_LIMITS.

**S1 = o-Terphenyl 40-140%**

S2 = 2-Fluorobiphenyl 40-140%

S3 = 1-Chlorooctadecane 40-140%

S4 = 2-Bromonaphthalene 40-140%

QC Limits (%)\* (Aqueous)

<u>LL to UL</u>	<u>40 to 140</u>	<u>40 to 140</u>	<u>40 to 140</u>	<u>40 to 140</u>
-----------------	------------------	------------------	------------------	------------------

QC Limits\* (Solid)

LL to UL           to        to        to        to

**Note:**

It is recommended that surrogate standard recoveries be monitored and documented on a continuing basis. At a minimum, when surrogate recovery from a sample, blank, or QC sample is less than 40% or more than 140%, check calculations to locate possible errors, check the fortifying standard solution for degradation, and check changes in instrument performance.

If the cause cannot be determined, reanalyze the sample unless one of the following exceptions applies:

- (1) Obvious interference is present on the chromatogram (e.g., unresolved complex mixture);
- (2) The surrogate exhibits high recovery and associated target analytes or hydrocarbon ranges are not detected in sample.

If a sample with a surrogate recovery outside of the acceptable range is not reanalyzed based on any of these aforementioned exceptions, this information must be noted on the data report form and discussed in the Executive Report. Analysis of the sample on dilution may diminish matrix-related surrogate recovery problems. This approach can be used as long as the reporting limits to evaluate applicable MCP standards can still be achieved with the dilution. If not, reanalysis without dilution must be performed.

## DATA REVIEW WORKSHEETS

All criteria were met   X    
Criteria were not met and/or see below           

### VII. A MATRIX SPIKE/MATRIX SPIKE DUPLICATE (MS/MSD)

This data is generated to determine long term precision and accuracy in the analytical method for various matrices. This data alone cannot be used to evaluate the precision and accuracy of individual samples.

At the request of the data user, and in consideration of sample matrices and data quality objectives, matrix spikes and matrix duplicates may be analyzed with every batch of 20 samples or less per matrix.

- **Matrix duplicate** - Matrix duplicates are prepared by analyzing one sample in duplicate. The purpose of the matrix duplicates is to determine the homogeneity of the sample matrix as well as analytical precision. The RPD of detected results in the matrix duplicate samples must not exceed 50 when the results are greater than 5x the reporting limit.
- The desired spiking level is 50% of the highest calibration standard. However, the total concentration in the MS (including the MS and native concentration in the unspiked sample) should not exceed 75% of the highest calibration standard in order for a proper evaluation to be performed. The purpose of the matrix spike is to determine whether the sample matrix contributes bias to the analytical results. The corrected concentrations of each analyte within the matrix spiking solution must be within 40 - 140% of the true value. Lower recoveries of n-nonane are permissible but must be noted in the narrative if <30%.

#### MS/MSD Recoveries and Precision Criteria

Sample ID: MC47204-2

Matrix/Level: Groundwater

List the %Rs, RPD of the compounds which do not meet the QC criteria.

MS OR MSD	COMPOUND	% R	RPD	QC LIMITS	ACTION

**Note:** MS/MSD analyzed with this data package. % recoveries and RPD within laboratory control limits.

## DATA REVIEW WORKSHEETS

All criteria were met   X    
 Criteria were not met and/or see below       

No action is taken on MS/MSD results alone to qualify the entire case. However, used informed professional judgment, the data reviewer may use the MS/MSD results in conjunction with other QC criteria and determine the need for some qualification of the data. In those instances where it can be determined that the results of the MS/MSD affect only the sample spiked, the qualification should be limited to this sample alone. However, it may be determined through the MS/MSD results that the laboratory is having a systematic problem in the analysis of one or more analytes, which affects the associated samples.

### 2. MS/MSD – Unspiked Compounds

List the concentrations of the unspiked compounds and determine the % RSDs of these compounds in the unspiked sample, matrix spike, and matrix spike duplicate.

COMPOUND	CONCENTRATION		MSD	%RPD	ACTION
	SAMPLE	MS			

Criteria: None specified, use %RSD  $\leq$  50 as professional judgment.

Actions:

If the % RSD > 50, qualify the results in the spiked sample as estimate (J).

If the % RSD is not calculable (NC) due to nondetect value in the sample, MS, and/or MSD, use professional judgment to qualify sample data.

A separate worksheet should be used for each MS/MSD pair.



## DATA REVIEW WORKSHEETS

All criteria were met   X    
Criteria were not met and/or see below

### VIII. LABORATORY CONTROL SAMPLE (LCS/LCSD) ANALYSIS

This data is generated to determine accuracy of the analytical method for various matrices.

## 1. LCS Recoveries Criteria

List the %R of compounds which do not meet the criteria

LCS ID	COMPOUND	% R	QC LIMIT	ACTION
--------	----------	-----	----------	--------

LCS RECOVERY WITHIN LABORATORY CONTROL LIMITS

**Criteria:**

- \* Refer to QAPP for specific criteria.
- \* The spike recovery must be between 40% and 140%. Lower recoveries of n-nonane are permissible. If the recovery of n-nonane is <30%, note the nonconformance in the executive narrative. RPD between LCS/LCSD must be < 25%.

**Actions:**

Actions on LCS recovery should be based on both the number of compounds that are outside the %R and RPD criteria and the magnitude of the exceedance of the criteria.

If the %R of the analyte is > UL, qualify all positive results (j) for the affected analyte in the associated samples and accept nondetects.

If the %R of the analyte is < LL, qualify all positive results (j) and reject (R) nondetects for the affected analyte in the associated samples.

If more than half the compounds in the LCS are not within the required recovery criteria, qualify all positive results as (J) and reject nondetects (R) for all target analyte(s) in the associated samples.

## 2. Frequency Criteria:

Where LCS analyzed at the required frequency and for each matrix (1 per 20 samples per matrix)? Yes or No.

If no, the data may be affected. Use professional judgment to determine the severity of the effect and qualify data accordingly. Discuss any actions below and list the samples affected. Discuss the actions below:

## DATA REVIEW WORKSHEETS

All criteria were met \_\_\_\_\_  
 Criteria were not met and/or see below N/A

### IX. FIELD/LABORATORY DUPLICATE PRECISION

Sample IDs: \_\_\_\_\_ - \_\_\_\_\_

Matrix: \_\_\_\_\_ - \_\_\_\_\_

Field/laboratory duplicate samples may be taken and analyzed as an indication of overall precision. These analyses measure both field and lab precision; therefore, the results may have more variability than laboratory duplicates which measures only laboratory performance. It is also expected that soil duplicate results will have a greater variance than water matrices due to difficulties associated with collecting identical field duplicate samples.

COMPOUND	SQL	SAMPLE CONC.	DUPLICATE CONC.	RPD	ACTION
No field/laboratory duplicate analyzed with this data package. MS/MSD % recoveries RPD used to assess precision. RPD within laboratory and validation guidance document criteria (+ 50 % RPD) for analytes concentration $\geq 5$ SQL.					

#### Criteria:

The project QAPP should be reviewed for project-specific information.  
 RPD  $\pm 30\%$  for aqueous samples, RPD  $\pm 50\%$  for solid samples if results are  $\geq$  SQL.  
 If both samples and duplicate are  $< 5$  SQL, the RPD criteria is doubled.

SQL = soil quantitation limit

#### Actions:

If both the sample and the duplicate results are nondetects (ND), the RPD is not calculable (NC). No action is needed.

Qualify as estimated positive results (J) and nondetects (UJ) for the compound that exceeded the above criteria.

If one sample result is not detected and the other is  $\geq 5$ x the SQL qualify (J/UJ).

**Note:** If SQLs for the sample and duplicate are significantly different, use professional judgment to determine if qualification is appropriate.

If one sample value is not detected and the other is  $< 5$ x the SQL, use professional judgment to determine if qualification is appropriate.

## DATA REVIEW WORKSHEETS

All criteria were met   X    
Criteria were not met and/or see below           

### XI. COMPOUND IDENTIFICATION

The compound identification evaluation is to verify that the laboratory correctly identified target analytes as well as tentatively identified compounds (TICs).

1. Verify that the target analytes were within the retention time windows.
  - Retention time windows must be re-established for each Target EPH Analyte each time a new GC column is installed, and must be verified and/or adjusted on a daily basis.
  - The n-nonane (n-C9) peak must be adequately resolved from the solvent front of the chromatographic run.
  - All surrogates must be adequately resolved from the Aliphatic Hydrocarbon and Aromatic Hydrocarbon standards.
  - For the purposes of this method, adequate resolution is assumed to be achieved if the height of the valley between two peaks is less than 25% of the average height of the two peaks.
  - The n-pentane (C5) and MtBE peaks must be adequately resolved from any solvent front that may be present on the FID and PID chromatograms, respectively.

- 1a. Aliphatic hydrocarbons range:
  - Determine the total area count for all peaks eluting 0.1 minutes before the retention time (Rt) for n-C9 and 0.01 minutes before the Rt for n-C19.
  - Determine the total area count for all peaks eluting 0.01 minutes before the Rt for n-C19 and 0.1 minutes after the Rt for n-C36.

Are the aliphatic hydrocarbons range properly determined? Yes? or No?

Comments:

- 1b. Aromatic hydrocarbons range:
  - Determine the total area count for all peaks eluting 0.1 minutes before the retention time (Rt) for naphthalene and 0.1 minutes after the Rt for benzo(g,h,i)perylene.
  - Determine the peak area count for the sample surrogate (OTP) and fractionation surrogate(s). Subtract these values from the collective area count value.

Are the aliphatic hydrocarbons range properly determined? Yes? or No?

Comments:

DATA REVIEW WORKSHEETS

All criteria were met   X    
Criteria were not met and/or see below           

2. If target analytes and/or TICs were not correctly identified, request that the laboratory resubmit the corrected data.
3. Breakthrough determination - Each sample (field and QC sample) must be evaluated for potential breakthrough on a sample specific basis by evaluating the % recovery of the fractionation surrogate (2-bromonaphthalene) and on a batch basis by quantifying naphthalene and 2-methylnaphthalene in both the aliphatic and aromatic fractions of the LCS and LCSD. **If either the concentration of naphthalene or 2-methylnaphthalene in the aliphatic fraction exceeds 5% of the total concentration for naphthalene or 2-methylnaphthalene in the LCS or LCSD, fractionation must be repeated on all archived batch extracts.**

**NOTE:** The total concentration of naphthalene or 2-methylnaphthalene in the LCS/LCSD pair includes the summation of the concentration detected in the aliphatic fraction and the concentration detected in the aromatic fraction.

  Comments: Concentration in the aliphatic fraction < 5% of the total  
  concentration for naphthalene and 2-methylnaphthalene

\_\_\_\_\_

\_\_\_\_\_

\_\_\_\_\_

4. **Fractionation Check Standard** – A fractionation check solution is prepared containing 14 alkanes and 17 PAHs at a nominal concentration of 200 ng/μl of each constituent. The Fractionation Check Solution must be used to evaluate the fractionation efficiency of each new lot of silica gel/cartridges, and establish the optimum hexane volume required to efficiently elute aliphatic hydrocarbons while not allowing significant aromatic hydrocarbon breakthrough. For each analyte contained in the fractionation check solution, excluding n-nonane, the Percent Recovery must be between 40 and 140%. A 30% Recovery is acceptable for n-nonane.

Is a fractionation check standard analyzed?

Yes? or No?

Comments: Not applicable.

## DATA REVIEW WORKSHEETS

All criteria were met   X    
Criteria were not met and/or see below           

### XII. QUANTITATION LIMITS AND SAMPLE RESULTS

The sample quantitation evaluation is to verify laboratory quantitation results.

In order to demonstrate the absence of aliphatic mass discrimination, the response ratio of C28 to C20 must be at least 0.85. If <0.85, this nonconformance must be noted in the laboratory case narrative.

The chromatograms of Continuing Calibration Standards for aromatics must be reviewed to ensure that there are no obvious signs of mass discrimination.

Is aliphatic mass discrimination observed in the sample? Yes? or No?

Is aromatic mass discrimination observed in the sample? Yes? or No?

1. In the space below, please show a minimum of one sample calculation:

JC47204-2MS          EPH (C11 – C22, Aromatics)          RF = 114,553

[ ] = (43529961)/(114,553)

[ ] = 380 ug/ml    Ok

JC47204-2MS          EPH (C19 – C36, Aliphatics)          RF = 72,594

[ ] = (15425606)/(72,594)

[ ] = 212.5 ug/ml    Ok

## DATA REVIEW WORKSHEETS

2. If requested, verify that the results were above the laboratory method detection limit (MDLs).
3. If dilutions performed, were the SQLs elevated accordingly by the laboratory? List the affected samples and dilution factor in the table below.

SAMPLE ID	DILUTION FACTOR	REASON FOR DILUTION

If dilution was not performed, estimate results (J) for the affected compounds. List the affected samples/compounds:

---

---