



# Analytical Laboratory

Analytical Lab  
Page 1 of 29

13339 Hagers Ferry Road  
Huntersville, NC 28078-7929  
McGuire Nuclear Complex - MG03A2  
Phone: 980-875-5245 Fax: 980-875-4349

## Order Summary Report

**Order Number:** J11080089

Project Name: WWTS - Biweekly

Customer Name(s): Bill Kennedy, Melonie Martin, Wayne Chapman, Tom Johnson

Customer Address: 3195 Pine Hall Rd  
Mailcode: Belews Steam Station  
Belews Creek, NC 28012

Lab Contact: Jason C Perkins Phone: 980-875-5348

**Report Authorized By:** \_\_\_\_\_ **Date:** 8/24/2011  
**(Signature)**

---

### Program Comments:

FGD BiMonthly - LL Hg

### Data Flags & Calculations:

Any analytical tests or individual analytes within a test flagged with a Qualifier indicate a deviation from the method quality system or quality control requirement. The qualifier description is found at the end of the Certificate of Analysis (sample results) under the qualifiers heading. All results are reported on a dry weight basis unless otherwise noted.

### Data Package:

This data package includes analytical results that are applicable only to the samples described in this narrative. An estimation of the uncertainty of measurement for the results in the report is available upon request. This report shall not be reproduced, except in full, without the written consent of the Analytical Laboratory. Please contact the Analytical laboratory with any questions. The order of individual sections within this report is as follows:

*Job Summary Report, Sample Identification, Technical Validation of Data Package, Analytical Laboratory Certificate of Analysis, Analytical Laboratory QC Reports, Sub-contracted Laboratory Results, Customer Specific Data Sheets, Reports & Documentation, Customer Database Entries, Test Case Narratives, Chain of Custody (COC)*

### Certification:

The Analytical Laboratory holds the following State Certifications : North Carolina (DENR) Certificate #248, South Carolina (DHEC) Laboratory ID # 99005. Contact the Analytical Laboratory for definitive information about the certification status of specific methods.

## Sample ID's & Descriptions:

Sample ID	Plant/Station	Collection Date and Time	Collected By	Sample Description
2011016915	BELEWS	10-Aug-11 8:00 AM	T. Johnson	FGD Purge Eff
2011016916	BELEWS	10-Aug-11 8:00 AM	T. Johnson	EQ TANK EFF.
2011016917	BELEWS	10-Aug-11 8:00 AM	T. Johnson	BIOREACTOR 1 INF.
2011016918	BELEWS	10-Aug-11 8:00 AM	T. Johnson	BIOREACTOR 2 INF.
2011016919	BELEWS	10-Aug-11 8:00 AM	T. Johnson	BIOREACTOR 2 EFF.
2011016920	BELEWS	02-Aug-11 2:15 PM	L.DAVIS	FILTER BLANK
2011016921	BELEWS	02-Aug-11 2:15 PM	L.DAVIS	Trip Blank
2011016934	BELEWS	10-Aug-11 12:20 PM	David Morris (Prism)	BIOREACTOR 1 INF.
2011016935	BELEWS	02-Aug-11 2:15 PM	L.DAVIS	HG BLANK BIOREACTOR 1 INF.
2011016936	BELEWS	10-Aug-11 12:30 PM	David Morris (Prism)	BIOREACTOR 2 INF.
2011016937	BELEWS	02-Aug-11 2:15 PM	L.DAVIS	Hg Blk BioReactor 2 Inf
2011016938	BELEWS	10-Aug-11 12:25 PM	David Morris (Prism)	BIOREACTOR 2 EFF.
2011016939	BELEWS	02-Aug-11 2:15 PM	L.DAVIS	Hg Blk BioReactor 2 Eff
13 Total Samples				

# Technical Validation Review

## Checklist:

COC and .pdf report are in agreement with sample totals and analyses (compliance programs and procedures).

☒ Yes

☐ No

All Results are less than the laboratory reporting limits.

☐ Yes

☒ No

All laboratory QA/QC requirements are acceptable.

☒ Yes

☐ No

The Vendor Laboratories have been qualified by the Analytical Laboratory

Yes

## Report Sections Included:

☒ Job Summary Report

☒ Sample Identification

☒ Technical Validation of Data Package

☒ Analytical Laboratory Certificate of Analysis

☐ Analytical Laboratory QC Report

☒ Sub-contracted Laboratory Results

☒ Customer Specific Data Sheets, Reports, & Documentation

☐ Customer Database Entries

☒ Chain of Custody

☐ Electronic Data Deliverable (EDD) Sent Separately

Reviewed By: Mary Ann Ogle

Date: 8/24/2011

**Certificate of Laboratory Analysis***This report shall not be reproduced, except in full.***Order # J11080089**

Site: FGD Purge Eff

Collection Date: 10-Aug-11 8:00 AM

**Sample #: 2011016915**

Matrix: OTHER

Analyte	Result	Units	Qualifiers	RDL	Method	Analysis Date/Time	Analyst
<b><u>INORGANIC IONS BY IC</u></b>							
Bromide	97	mg/L		10	EPA 300.0	11-Aug-11 20:40	CLEEMAN
<b><u>MERCURY (COLD VAPOR) IN WATER</u></b>							
Mercury (Hg)	244	ug/L		5	EPA 245.1	12-Aug-11 09:37	AGIBBS
<b><u>TOTAL RECOVERABLE METALS BY ICP</u></b>							
Boron (B)	168	mg/L		0.5	EPA 200.7	15-Aug-11 12:28	DJSULL1
<b><u>DISSOLVED METALS BY ICP-MS</u></b>							
Selenium (Se)	220	ug/L		10	EPA 200.8	15-Aug-11 12:11	KRICHAR
<b><u>TOTAL RECOVERABLE METALS BY ICP-MS</u></b>							
Arsenic (As)	303	ug/L		10	EPA 200.8	17-Aug-11 10:46	MHH7131
Chromium (Cr)	329	ug/L		10	EPA 200.8	17-Aug-11 10:46	MHH7131
Copper (Cu)	259	ug/L		10	EPA 200.8	17-Aug-11 10:46	MHH7131
Nickel (Ni)	260	ug/L		10	EPA 200.8	17-Aug-11 10:46	MHH7131
Selenium (Se)	7470	ug/L		20	EPA 200.8	17-Aug-11 10:46	MHH7131
Silver (Ag)	< 10	ug/L		10	EPA 200.8	17-Aug-11 10:46	MHH7131
Zinc (Zn)	466	ug/L		20	EPA 200.8	17-Aug-11 10:46	MHH7131
<b><u>SELENIUM SPECIATION</u></b>							
Vendor Parameter	Complete				V_AS&C		
<b><u>TOTAL DISSOLVED SOLIDS</u></b>							
TDS	14000	mg/L		10	SM2540C	12-Aug-11 14:30	TJA7067

Site: EQ TANK EFF.

Collection Date: 10-Aug-11 8:00 AM

**Sample #: 2011016916**

Matrix: OTHER

Analyte	Result	Units	Qualifiers	RDL	Method	Analysis Date/Time	Analyst
<b><u>MERCURY (COLD VAPOR) IN WATER</u></b>							
Mercury (Hg)	317	ug/L		5	EPA 245.1	12-Aug-11 09:45	AGIBBS
<b><u>TOTAL RECOVERABLE METALS BY ICP</u></b>							
Boron (B)	179	mg/L		0.5	EPA 200.7	15-Aug-11 12:31	DJSULL1
<b><u>DISSOLVED METALS BY ICP-MS</u></b>							
Selenium (Se)	326	ug/L		10	EPA 200.8	15-Aug-11 12:14	KRICHAR
<b><u>TOTAL RECOVERABLE METALS BY ICP-MS</u></b>							
Arsenic (As)	345	ug/L		10	EPA 200.8	17-Aug-11 10:46	MHH7131

# Certificate of Laboratory Analysis

Analytical Lab  
Page 5 of 29

*This report shall not be reproduced, except in full.*

Order # J11080089

Site: EQ TANK EFF.

Collection Date: 10-Aug-11 8:00 AM

Sample #: 2011016916

Matrix: OTHER

Analyte	Result	Units	Qualifiers	RDL	Method	Analysis Date/Time	Analyst
<b><u>TOTAL RECOVERABLE METALS BY ICP-MS</u></b>							
Chromium (Cr)	412	ug/L		10	EPA 200.8	17-Aug-11 10:46	MHH7131
Copper (Cu)	315	ug/L		10	EPA 200.8	17-Aug-11 10:46	MHH7131
Nickel (Ni)	295	ug/L		10	EPA 200.8	17-Aug-11 10:46	MHH7131
Selenium (Se)	9530	ug/L		20	EPA 200.8	17-Aug-11 10:46	MHH7131
Silver (Ag)	< 10	ug/L		10	EPA 200.8	17-Aug-11 10:46	MHH7131
Zinc (Zn)	541	ug/L		20	EPA 200.8	17-Aug-11 10:46	MHH7131

Site: BIOREACTOR 1 INF.

Collection Date: 10-Aug-11 8:00 AM

Sample #: 2011016917

Matrix: OTHER

Analyte	Result	Units	Qualifiers	RDL	Method	Analysis Date/Time	Analyst
<b><u>TOTAL RECOVERABLE METALS BY ICP</u></b>							
Boron (B)	167	mg/L		0.5	EPA 200.7	15-Aug-11 12:35	DJSULL1
<b><u>DISSOLVED METALS BY ICP-MS</u></b>							
Selenium (Se)	210	ug/L		5	EPA 200.8	15-Aug-11 12:17	KRICAR
<b><u>TOTAL RECOVERABLE METALS BY ICP-MS</u></b>							
Arsenic (As)	< 10	ug/L		10	EPA 200.8	17-Aug-11 10:43	MHH7131
Chromium (Cr)	< 10	ug/L		10	EPA 200.8	17-Aug-11 10:43	MHH7131
Copper (Cu)	< 10	ug/L		10	EPA 200.8	17-Aug-11 10:43	MHH7131
Nickel (Ni)	< 10	ug/L		10	EPA 200.8	17-Aug-11 10:43	MHH7131
Selenium (Se)	242	ug/L		10	EPA 200.8	17-Aug-11 10:43	MHH7131
Silver (Ag)	< 10	ug/L		10	EPA 200.8	17-Aug-11 10:43	MHH7131
Zinc (Zn)	< 20	ug/L		20	EPA 200.8	17-Aug-11 10:43	MHH7131

## **SELENIUM SPECIATION**

Vendor Parameter Complete V\_AS&C

Site: BIOREACTOR 2 INF.

Collection Date: 10-Aug-11 8:00 AM

Sample #: 2011016918

Matrix: OTHER

Analyte	Result	Units	Qualifiers	RDL	Method	Analysis Date/Time	Analyst
<b><u>TOTAL RECOVERABLE METALS BY ICP</u></b>							
Boron (B)	168	mg/L		0.5	EPA 200.7	15-Aug-11 12:39	DJSULL1
<b><u>TOTAL RECOVERABLE METALS BY ICP-MS</u></b>							
Arsenic (As)	< 10	ug/L		10	EPA 200.8	17-Aug-11 10:46	MHH7131
Chromium (Cr)	< 10	ug/L		10	EPA 200.8	17-Aug-11 10:46	MHH7131
Copper (Cu)	< 10	ug/L		10	EPA 200.8	17-Aug-11 10:46	MHH7131
Nickel (Ni)	< 10	ug/L		10	EPA 200.8	17-Aug-11 10:46	MHH7131

# Certificate of Laboratory Analysis

Analytical Lab  
Page 6 of 29

*This report shall not be reproduced, except in full.*

Order # J11080089

Site: BIOREACTOR 2 INF.

Collection Date: 10-Aug-11 8:00 AM

Sample #: 2011016918

Matrix: OTHER

Analyte	Result	Units	Qualifiers	RDL	Method	Analysis Date/Time	Analyst
<b>TOTAL RECOVERABLE METALS BY ICP-MS</b>							
Selenium (Se)	52.1	ug/L		10	EPA 200.8	17-Aug-11 10:46	MHH7131
Silver (Ag)	< 10	ug/L		10	EPA 200.8	17-Aug-11 10:46	MHH7131
Zinc (Zn)	< 20	ug/L		20	EPA 200.8	17-Aug-11 10:46	MHH7131

Site: BIOREACTOR 2 EFF.

Collection Date: 10-Aug-11 8:00 AM

Sample #: 2011016919

Matrix: OTHER

Analyte	Result	Units	Qualifiers	RDL	Method	Analysis Date/Time	Analyst
<b>INORGANIC IONS BY IC</b>							
Bromide	95	mg/L		10	EPA 300.0	11-Aug-11 20:55	CLEEMAN
<b>MERCURY (COLD VAPOR) IN WATER</b>							
Mercury (Hg)	< 1	ug/L		1	EPA 245.1	12-Aug-11 09:47	AGIBBS
<b>TOTAL RECOVERABLE METALS BY ICP</b>							
Boron (B)	173	mg/L		0.5	EPA 200.7	15-Aug-11 12:43	DJSULL1
<b>TOTAL RECOVERABLE METALS BY ICP-MS</b>							
Arsenic (As)	< 5	ug/L		5	EPA 200.8	17-Aug-11 10:46	MHH7131
Chromium (Cr)	< 5	ug/L		5	EPA 200.8	17-Aug-11 10:46	MHH7131
Copper (Cu)	< 5	ug/L		5	EPA 200.8	17-Aug-11 10:46	MHH7131
Nickel (Ni)	< 5	ug/L		5	EPA 200.8	17-Aug-11 10:46	MHH7131
Selenium (Se)	5.01	ug/L		5	EPA 200.8	17-Aug-11 10:46	MHH7131
Silver (Ag)	< 5	ug/L		5	EPA 200.8	17-Aug-11 10:46	MHH7131
Zinc (Zn)	< 10	ug/L		10	EPA 200.8	17-Aug-11 10:46	MHH7131

## SELENIUM SPECIATION

Vendor Parameter Complete V\_AS&C

Site: FILTER BLANK

Collection Date: 02-Aug-11 2:15 PM

Sample #: 2011016920

Matrix: OTHER

Analyte	Result	Units	Qualifiers	RDL	Method	Analysis Date/Time	Analyst
<b>DISSOLVED METALS BY ICP-MS</b>							
Selenium (Se)	< 1	ug/L		1	EPA 200.8	15-Aug-11 11:51	KRICHR

Site: Trip Blank

Collection Date: 02-Aug-11 2:15 PM

Sample #: 2011016921

Matrix: OTHER

Analyte	Result	Units	Qualifiers	RDL	Method	Analysis Date/Time	Analyst
---------	--------	-------	------------	-----	--------	--------------------	---------

# Certificate of Laboratory Analysis

Analytical Lab  
Page 7 of 29

*This report shall not be reproduced, except in full.*

Order # J11080089

Site: Trip Blank

Collection Date: 02-Aug-11 2:15 PM

Sample #: 2011016921

Matrix: OTHER

Analyte	Result	Units	Qualifiers	RDL	Method	Analysis Date/Time	Analyst
<b>TOTAL RECOVERABLE METALS BY ICP</b>							
Boron (B)	< 0.05	mg/L		0.05	EPA 200.7	15-Aug-11 12:24	DJSULL1
<b>TOTAL RECOVERABLE METALS BY ICP-MS</b>							
Arsenic (As)	< 1	ug/L		1	EPA 200.8	17-Aug-11 10:27	MHH7131
Chromium (Cr)	< 1	ug/L		1	EPA 200.8	17-Aug-11 10:27	MHH7131
Copper (Cu)	< 1	ug/L		1	EPA 200.8	17-Aug-11 10:27	MHH7131
Nickel (Ni)	< 1	ug/L		1	EPA 200.8	17-Aug-11 10:27	MHH7131
Selenium (Se)	< 1	ug/L		1	EPA 200.8	17-Aug-11 10:27	MHH7131
Silver (Ag)	< 1	ug/L		1	EPA 200.8	17-Aug-11 10:27	MHH7131
Zinc (Zn)	< 2	ug/L		2	EPA 200.8	17-Aug-11 10:27	MHH7131

## SELENIUM SPECIATION

Vendor Parameter Complete V\_AS&C

Site: BIOREACTOR 1 INF.

Collection Date: 10-Aug-11 12:20 PM

Sample #: 2011016934

Matrix: OTHER

Analyte	Result	Units	Qualifiers	RDL	Method	Analysis Date/Time	Analyst
<b>MERCURY 1631</b>							
Vendor Parameter	Complete				V_BRAND		

Site: HG BLANK BIOREACTOR 1 INF.

Collection Date: 02-Aug-11 2:15 PM

Sample #: 2011016935

Matrix: OTHER

Analyte	Result	Units	Qualifiers	RDL	Method	Analysis Date/Time	Analyst
<b>MERCURY 1631</b>							
Vendor Parameter	Complete				V_BRAND		

Site: BIOREACTOR 2 INF.

Collection Date: 10-Aug-11 12:30 PM

Sample #: 2011016936

Matrix: OTHER

Analyte	Result	Units	Qualifiers	RDL	Method	Analysis Date/Time	Analyst
<b>MERCURY 1631</b>							
Vendor Parameter	Complete				V_BRAND		

Site: Hg Blk BioReactor 2 Inf

Collection Date: 02-Aug-11 2:15 PM

Sample #: 2011016937

Matrix: OTHER

Analyte	Result	Units	Qualifiers	RDL	Method	Analysis Date/Time	Analyst
---------	--------	-------	------------	-----	--------	--------------------	---------

# Certificate of Laboratory Analysis

Analytical Lab  
Page 8 of 29

*This report shall not be reproduced, except in full.*

**Order # J11080089**

Site: Hg Blk BioReactor 2 Inf

Collection Date: 02-Aug-11 2:15 PM

**Sample #:** 2011016937

Matrix: OTHER

Analyte	Result	Units	Qualifiers	RDL	Method	Analysis Date/Time	Analyst
<b><u>MERCURY 1631</u></b>							
Vendor Parameter	Complete				V_BRAND		

Site: BIOREACTOR 2 EFF.

Collection Date: 10-Aug-11 12:25 PM

**Sample #:** 2011016938

Matrix: OTHER

Analyte	Result	Units	Qualifiers	RDL	Method	Analysis Date/Time	Analyst
<b><u>MERCURY 1631</u></b>							
Vendor Parameter	Complete				V_BRAND		

Site: Hg Blk BioReactor 2 Eff

Collection Date: 02-Aug-11 2:15 PM

**Sample #:** 2011016939

Matrix: OTHER

Analyte	Result	Units	Qualifiers	RDL	Method	Analysis Date/Time	Analyst
<b><u>MERCURY 1631</u></b>							
Vendor Parameter	Complete				V_BRAND		





**APPLIED SPECIATION  
AND CONSULTING, LLC**

18804 Northcreek Parkway Bothell, WA, 98011  
Tel: (425) 483-3300 Fax: (425) 483-9818  
[www.appliedspeciation.com](http://www.appliedspeciation.com)

August 22, 2011

Jay Perkins  
Duke Energy Analytical Laboratory  
Mail Code MGO3A2 (Building 7405)  
13339 Hagers Ferry Rd.  
Huntersville, NC 28078  
(704) 875-5245

Project: Belews – FGD WWTS (Bi-Monthly Sampling) (LIMS #J11080089)

Dear Mr. Perkins,

Attached is the report associated with four (4) aqueous samples submitted for selenium speciation analysis on August 15, 2011. The samples were received in a sealed cooler at  $-0.4^{\circ}\text{C}$  on August 16, 2011. Selenium speciation analysis was performed via ion chromatography inductively coupled plasma dynamic reaction cell mass spectrometry (IC-ICP-DRC-MS). Any issues associated with the analysis are addressed in the following report.

If you have any questions, please feel free to contact me at your convenience.

Sincerely,

A handwritten signature in black ink that reads "Ben Wozniak".

Ben Wozniak  
Project Manager  
Applied Speciation and Consulting, LLC

Applied Speciation and Consulting, LLC

Report prepared for:

Jay Perkins  
Duke Energy Analytical Laboratory  
Mail Code MGO3A2 (Building 7405)  
13339 Hagers Ferry Rd.  
Huntersville, NC 28078

Project: Belews – FGD WWTS (Bi-Monthly Sampling) (LIMS #J11080089)

August 22, 2011

## 1. Sample Reception

Four (4) aqueous samples in 125mL HDPE bottles (provided by Applied Speciation and Consulting) were submitted for selenium speciation analysis on August 15, 2011. The samples were received on August 16, 2011 in a sealed container at -0.4°C.

The samples were received in a laminar flow clean hood void of trace metals contamination and ultra-violet radiation. Upon reception, the samples were designated discrete sample identifiers. An aliquot of each sample was filtered (0.45µm) and these filtrates were stored in a secure, monitored cryofreezer (maintained at a temperature of -80°C) until selenium speciation analysis could be performed via ion chromatography inductively coupled plasma dynamic reaction cell mass spectrometry (IC-ICP-DRC-MS).

## 2. Sample Preparation

All sample preparation is performed in laminar flow clean hoods known to be free from trace metals contamination. All applied water for dilutions and sample preservatives are monitored for contamination to account for any biases associated with the sample results.

Selenium Speciation Analysis by IC-ICP-DRC-MS Prior to analysis, an aliquot of each sample was filtered with a syringe filter (0.45µm) and injected directly into a sealed autosampler vial. No further sample preparation was performed as any chemical alteration of the samples may shift the equilibrium of the system resulting in changes in speciation ratios.

## 3. Sample Analysis

All sample analysis is preceded by a minimum of a five-point calibration curve spanning the entire concentration range of interest. Calibration curves are performed at the beginning of

each analytical day. All calibration curves, associated with each species of interest, are standardized by linear regression resulting in a response factor. All sample results are **instrument blank corrected** to account for any operational biases associated with the analytical platform.

Prior to sample analysis, all calibration curves are verified using second source standards which are identified as initial calibration verification standards (ICV).

Ongoing instrument performance is identified by the analysis of continuing calibration verification standards (CCV) and continuing calibration blanks (CCB) at a minimal interval of every ten analytical runs.

*Selenium Speciation Analysis by IC-ICP-DRC-MS* All samples for selenium speciation analysis were analyzed by ion chromatography inductively coupled plasma dynamic reaction cell mass spectrometry (IC-ICP-DRC-MS) on August 17, 2011. An aliquot of each sample is injected onto an anion exchange column and mobilized by a basic ( $\text{pH} > 7$ ) gradient. The eluting selenium species are then introduced into a radio frequency (RF) plasma where energy-transfer processes cause desolvation, atomization, and ionization. The ions are extracted from the plasma through a differentially-pumped vacuum interface and travel through a pressurized chamber (DRC) containing a specific reactive gas which preferentially reacts with interfering ions of the same target mass to charge ratios ( $m/z$ ). A solid-state detector detects ions transmitted through the mass analyzer and the resulting current is processed by a data handling system.

Retention times for each eluting species are compared to known standards for species identification.

#### 4. Analytical Issues


The overall analyses went well and no significant analytical issues were encountered. All quality control parameters associated with these samples were within acceptance limits.

The estimated method detection limits (eMDLs) for selenite, selenate, and selenocyanate are generated from replicate analyses of the lowest standard in the calibration curve. Not all selenium species are present in preparation blanks; therefore, eMDL calculations based on preparation blanks are artificially biased low.

The eMDL for methylseleninic acid and selenomethionine is calculated from the average eMDL of selenite, selenate, and selenocyanate. The calibration does not contain methylseleninic acid or selenomethionine due to impurities in these standards which would bias the results for other selenium species.

If you have any questions or concerns regarding this report, please feel free to contact me.

Sincerely,

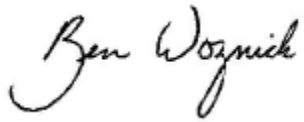A handwritten signature in black ink that reads "Ben Wozniak". The signature is written in a cursive style with a large, looping 'B' and a trailing flourish.

Ben Wozniak  
Project Manager  
Applied Speciation and Consulting, LLC

Selenium Speciation Results for Duke Energy  
Project Name: Belews – FGD WWTS (Bi-Monthly Sampling)  
Contact: Jay Perkins  
LIMS #J11080089

Date: August 22, 2011  
Report Generated by: Ben Wozniak  
Applied Speciation and Consulting, LLC

**Sample Results**

Sample ID	Se(IV)	Se(VI)	SeCN	MeSe(IV)	SeMe	Unknown Se Species (n)
FGD Purge Eff	203	138	ND (<3.9)	ND (<4.6)	ND (<4.6)	0 (0)
BioReactor 1 Inf	99.6	122	3.76	2.7	ND (<1.1)	5.4 (2)
BioReactor 2 Eff	ND (<1.6)	ND (<0.90)	ND (<0.98)	ND (<1.1)	ND (<1.1)	0 (0)
Metals Trip Blk	ND (<0.31)	ND (<0.18)	ND (<0.20)	ND (<0.23)	ND (<0.23)	0 (0)

All results reflect the applied dilution and are reported in µg/L

ND = Not detected at the applied dilution

SeCN = Selenocyanate

MeSe(IV) = Methylseleninic acid

SeMe = Selenomethionine

Unknown Se Species = Total concentration of all unknown Se species observed by IC-ICP-MS

n = number of unknown Se species observed

Selenium Speciation Results for Duke Energy  
Project Name: Belews – FGD WWTS (Bi-Monthly Sampling)  
Contact: Jay Perkins  
LIMS #J11080089

Date: August 22, 2011  
Report Generated by: Ben Wozniak  
Applied Speciation and Consulting, LLC

**Quality Control Summary - Preparation Blank Summary**

Analyte (µg/L)	PBW1	PBW2	PBW3	PBW4	Mean	StdDev	eMDL *	eMDL 10x	eMDL 50x	eMDL 200x
Se(IV)	0.000	0.000	0.000	0.000	0.000	0.000	0.031	0.31	1.6	6.2
Se(VI)	0.000	0.000	0.000	0.000	0.000	0.000	0.018	0.18	0.90	3.6
SeCN	0.000	0.000	0.000	0.000	0.000	0.000	0.020	0.20	0.98	3.9
MeSe(IV)	0.000	0.000	0.000	0.000	0.000	0.000	0.023	0.23	1.1	4.6
SeMe	0.000	0.000	0.000	0.000	0.000	0.000	0.023	0.23	1.1	4.6

eMDL = Estimated Method Detection Limit

\*Please see narrative regarding eMDL calculations

**Quality Control Summary - Certified Reference Materials**

Analyte (µg/L)	CRM	True Value	Result	Recovery
Se(IV)	LCS	9.57	9.759	102.0
Se(VI)	LCS	9.48	9.480	100.0
SeCN	LCS	8.92	8.929	100.1
MeSe(IV)	LCS	6.47	6.625	102.4
SeMe	LCS	9.32	9.408	100.9

Selenium Speciation Results for Duke Energy  
Project Name: Belews – FGD WWTS (Bi-Monthly Sampling)  
Contact: Jay Perkins  
LIMS #J11080089

Date: August 22, 2011  
Report Generated by: Ben Wozniak  
Applied Speciation and Consulting, LLC

**Quality Control Summary - Matrix Duplicates**

Analyte (µg/L)	Sample ID	Rep 1	Rep 2	Mean	RPD
Se(IV)	BioReactor 2 Eff	ND (<1.6)	ND (<1.6)	NC	NC
Se(VI)	BioReactor 2 Eff	ND (<0.90)	ND (<0.90)	NC	NC
SeCN	BioReactor 2 Eff	ND (<0.98)	ND (<0.98)	NC	NC
MeSe(IV)	BioReactor 2 Eff	ND (<1.1)	ND (<1.1)	NC	NC
SeMe	BioReactor 2 Eff	ND (<1.1)	ND (<1.1)	NC	NC

ND = Not detected at the applied dilution

NC = Value was not calculated due to one or more concentrations below the eMDL

**Quality Control Summary - Matrix Spike/ Matrix Spike Duplicate**

Analyte (µg/L)	Sample ID	Spike Conc	MS Result	Recovery	Spike Conc	MSD Result	Recovery	RPD
Se(IV)	BioReactor 2 Eff	278.0	275.8	99.2	278.0	274.2	98.6	0.6
Se(VI)	BioReactor 2 Eff	252.3	264.5	104.9	252.3	265.8	105.4	0.5
SeCN	BioReactor 2 Eff	228.8	223.7	97.8	228.8	225.1	98.4	0.6



# CHAIN OF CUSTODY RECORD AND ANALYSIS REQUEST FORM

Analytical Lab  
Page 18 of 29  
Page 1 of 2  
DISTRIBUTION  
ORIGINAL to LAB,  
COPY to CLIENT



## Duke Energy Analytical Laboratory

Mail Code MGO3A2 (Building 7405)  
13339 Hagers Ferry Rd  
Huntersville, N. C. 28078  
(704) 875-5245  
Fax: (704) 875-4349

## Analytical Laboratory Use Only

ORDER# <b>11080089</b>	MATRIX: OTHER	Samples Originating From	NC SC
Logged By <b>Cpk</b>	Date & Time <b>8-10-11 1546</b>	SAMPLE PROGRAM Water _____ Ground NPDES Drinking Water UST _____ RCRA Waste _____	
Ve <b>AS&amp;C</b>	PO#133241	Cooler Temp (C) <b>&lt; 1</b>	
MR #		15 Preserv.: 1=HCL 2=H <sub>2</sub> SO <sub>4</sub> 3=HNO <sub>3</sub> 4=Ice 5=None	

1) Project Name <b>Belews - FGD</b>	2) Phone No:
<b>WWTS Bi-Monthly Sampling)</b>	
2) Client: <b>Bill Kennedy, Melonie Martin, Wayne Chapman, Tom Johnson **</b>	4) Fax No:
5) Business Unit:	6) Process:
8) Oper. Unit:	9) Res. Type:
	10) Reso. Center:
	Mail Code:

Customer to complete all  
appropriate non-shaded areas.

Sampling conducted: 2nd and 4th Wednesday

Se Speciation Bottle ID	13 Sample Description or ID	Date	Time	Signature	17 Comp.	18 Grab	TDS	Hg - 245.1	Br (Dionex)	Metals*	Se, soluble (no dig.)	Se, speciation - vendor to AS&C (Important to place filled bottle back into both baggies)
B10871	FGD Purge Eff	8/10	0800	Tom Johnson		✓	1	1	1	1		1
	EQ Tank Eff.	8/10	0800	Tom Johnson		✓		1	1	1		1
B10859	BioReactor 1 Inf	8/10	0800		✓					1	1	1
	BioReactor 2 Inf	8/10	0800			✓				1		
B11239	BioReactor 2 Eff	8/10	0800		✓			1	1	1		1
	Filter Blk	8/2/11	1415	R. Davis		✓				1		1
B11716	Metals Trip Blk	8/2/11	1415	R. Davis		✓				1		1

Filtering of the Se is performed in the field please provide a filter blank too.

Customer to sign & date below - fill out from left to right.

1) Relinquished By <b>Tom Johnson</b>	Date/Time <b>8/10/11 0900</b>	2) Accepted By <b>Tom Johnson</b>	Date/Time <b>8-10-11 1215</b>
3) Relinquished By <b>David Mon</b>	Date/Time <b>8-10-11 1525</b>	4) Accepted By <b>Candy Knox</b>	Date/Time <b>8-10-11 1525</b>
5) Relinquished By <b>Cpk</b>	Date/Time <b>8-11-11 1300</b>	6) Accepted By <b>Nancy Cullinan</b>	Date/Time <b>8/10/11 9:15</b>
7) Relinquished By <b>Cpk</b>	Date/Time <b>8-11-11</b>	8) Accepted By <b>Nancy Cullinan</b>	Date/Time <b>8/10/11 9:15</b>
9) Seal/Locked By <b>Cpk</b>	Date/Time <b>8-11-11</b>	10) Seal/Lock Opened By	Date/Time
11) Seal/Locked By	Date/Time	12) Seal/Lock Opened By	Date/Time
Comments			

Customer, IMPORTANT!  
Please indicate desired turnaround.

## 22 Requested Turnaround

14 Days \_\_\_\_\_

\* 7 Days \_\_\_\_\_

\* 48 Hr \_\_\_\_\_

\* Other \_\_\_\_\_

\* Add. Cost Will Apply

\* B by ICP

As, Cr, Cu, Ni, Se, Ag, Zn by IMS

Digestions = TRM

thomas.d.johnson@siemens.com



August 23, 2011

Duke Energy  
ATTN: Jay Perkins  
Scientific Support-Laboratory  
13339 Hagers Ferry Road  
Huntersville NC 28078  
jcperkins@duke-energy.com  
labcustomer@duke-energy.com

RE: Project DUK-HV1101

Client Project: J11080089

Dear Mr. Perkins,


On August 11, 2011, Brooks Rand Labs (BRL) received three (3) flue gas desulfurization (FGD) waste water samples and three (3) corresponding field blanks. Samples were logged-in for total mercury (Hg) analysis. All samples were received, prepared, analyzed, and stored according to BRL SOPs and EPA methodology.

The results were blank-corrected as described in the calculations section of the applicable SOP(s) and may be evaluated using adjusted reporting limits to account for sample aliquot size. Please refer to the *Sample Results* page for sample-specific detection limits and other details. All other quality assurance criteria were satisfied and, aside from concentration qualifiers, all data was reported without additional qualification.

BRL, an accredited laboratory, certifies the reported results of all analyses for which BRL is NELAP accredited meet all NELAP requirements. For more details, see the *Report Information* page of the report.

Please feel free to contact me if you have any questions regarding this report.

Sincerely,

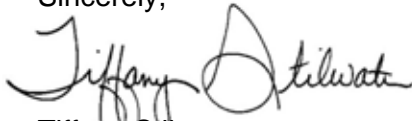

Tiffany Stilwater  
Project Manager  
tiffany@brooksrn.com

## Report Information

### Laboratory Accreditation

BRL is accredited by the *National Environmental Laboratory Accreditation Program* (NELAP) through the State of Florida Department of Health, Bureau of Laboratories (E87982) and is certified to perform many environmental analyses. BRL is also certified by many other states to perform environmental analyses. For a current list of our accreditations/certifications, please visit our website at <http://www.brooksrand.com/default.asp?contentID=586>. Results reported relate only to the samples listed in the report.

### Field Quality Control Samples

Please be notified that certain EPA methods require the collection of field quality control samples of an appropriate type and frequency; failure to do so is considered a deviation from some methods and for compliance purposes should only be done with the approval of regulatory authorities. Please see the specific EPA methods for details regarding required field quality control samples.

### Common Abbreviations


<b>BLK</b>	method blank	<b>MS</b>	matrix spike
<b>BRL</b>	Brooks Rand Labs	<b>MSD</b>	matrix spike duplicate
<b>BS</b>	laboratory fortified blank	<b>ND</b>	non-detect
<b>CAL</b>	calibration standard	<b>NR</b>	non-reportable
<b>CCV</b>	continuing calibration verification	<b>PS</b>	post preparation spike
<b>COC</b>	chain of custody record	<b>REC</b>	percent recovery
<b>CRM</b>	certified reference material	<b>RPD</b>	relative percent difference
<b>D</b>	dissolved fraction	<b>RSD</b>	relative standard deviation
<b>DUP</b>	duplicate	<b>SCV</b>	secondary calibration verification
<b>ICV</b>	initial calibration verification	<b>SOP</b>	standard operating procedure
<b>MDL</b>	method detection limit	<b>SRM</b>	standard reference material
<b>MRL</b>	method reporting limit	<b>T</b>	total recoverable fraction

### Definition of Data Qualifiers

(Effective 9/23/09)

<b>B</b>	Detected by the instrument, the result is > the MDL but ≤ the MRL. Result is reported and considered an estimate.
<b>E</b>	An estimated value due to the presence of interferences. A full explanation is presented in the narrative.
<b>H</b>	Holding time and/or preservation requirements not met. Result is estimated.
<b>J</b>	Estimated value. A full explanation is presented in the narrative.
<b>J-M</b>	Duplicate precision (RPD) for associated QC sample was not within acceptance criteria. Result is estimated.
<b>J-N</b>	Spike recovery for associated QC sample was not within acceptance criteria. Result is estimated.
<b>M</b>	Duplicate precision (RPD) was not within acceptance criteria. Result is estimated.
<b>N</b>	Spike recovery was not within acceptance criteria. Result is estimated.
<b>R</b>	Rejected, unusable value. A full explanation is presented in the narrative.
<b>U</b>	Result is ≤ the MDL or client requested reporting limit (CRRL). Result reported as the MDL or CRRL.
<b>X</b>	Result is not BLK-corrected and is within 10x the absolute value of the highest detectable BLK in the batch. Result is estimated.

These qualifiers are based on those previously utilized by Brooks Rand, Ltd., those found in the EPA SOW ILM03.0, Exhibit B, Section III, pg. B-18, and the USEPA Laboratory Data Validation Functional Guidelines for Evaluating Inorganic Analyses; USEPA; July 2002. These supersede all previous qualifiers ever employed by BRL.

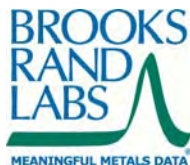


## Sample Information

Sample	Lab ID	Report Matrix	Type	Sampled	Received
BioReactor 1 Inf	1134009-01	FGD Wastewater	Sample	08/10/2011	08/16/2011
Hg Blk BioReactor 1 Inf	1134009-02	DIW	Field Blank	08/02/2011	08/16/2011
BioReactor 2 Inf	1134009-03	FGD Wastewater	Sample	08/10/2011	08/16/2011
Hg Blk BioReactor 2 Inf	1134009-04	DIW	Field Blank	08/02/2011	08/16/2011
BioReactor 2 Eff	1134009-05	FGD Wastewater	QC Sample	08/10/2011	08/16/2011
Hg Blk BioReactor 2 Eff	1134009-06	DIW	Field Blank	08/02/2011	08/16/2011

## Batch Summary

Analyte	Lab Matrix	Method	Prepared	Analyzed	Batch	Sequence
Hg	Water	EPA 1631	08/18/2011	08/22/2011	B111247	1100569

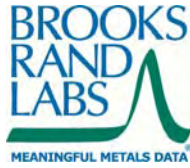

## Sample Results

Sample	Analyte	Report Matrix	Fraction	Result	Qualifier	MDL	MRL	Unit	Batch	Sequence
<b>BioReactor 1 Inf</b>										
1134009-01	Hg	FGD Wastewater	T	60.2		3.03	8.08	ng/L	B111247	1100569
<b>BioReactor 2 Eff</b>										
1134009-05	Hg	FGD Wastewater	T	31.1		0.78	2.07	ng/L	B111247	1100569
<b>BioReactor 2 Inf</b>										
1134009-03	Hg	FGD Wastewater	T	49.9		3.03	8.08	ng/L	B111247	1100569
<b>Hg Blk BioReactor 1 Inf</b>										
1134009-02	Hg	DIW	T	0.15	U	0.15	0.40	ng/L	B111247	1100569
<b>Hg Blk BioReactor 2 Eff</b>										
1134009-06	Hg	DIW	T	0.15	U	0.15	0.41	ng/L	B111247	1100569
<b>Hg Blk BioReactor 2 Inf</b>										
1134009-04	Hg	DIW	T	0.15	U	0.15	0.41	ng/L	B111247	1100569

## Accuracy & Precision Summary

Batch: B111247  
Lab Matrix: Water  
Method: EPA 1631


Sample	Analyte	Native	Spike	Result	Units	REC & Limits	RPD & Limits
B111247-SRM1	Certified Reference Material (1133002, NIST 1641d 1000x dilution)						
	Hg		15.68	14.68	ng/L	94% 85-115	
B111247-MS1	Matrix Spike (1134009-05)						
	Hg	31.14	157.1	214.1	ng/L	116% 71-125	
B111247-MSD1	Matrix Spike Duplicate (1134009-05)						
	Hg	31.14	153.3	204.7	ng/L	113% 71-125	5% 24

## Method Blanks & Reporting Limits

Batch: B111247  
Matrix: Water  
Method: EPA 1631  
Analyte: Hg

Sample	Result	Units
B111247-BLK1	0.10	ng/L
B111247-BLK2	0.05	ng/L
B111247-BLK3	0.08	ng/L
B111247-BLK4	0.06	ng/L
Average: 0.07		Standard Deviation: 0.02
Limit: 0.50		Limit: 0.10
		MDL: 0.15
		MRL: 0.41

**Project ID:** DUK-HV1101  
**PM:** Tiffany Stilwater

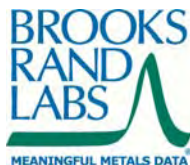

Analytical Lab  
Page 22 of 29  
**Client PM:** Jay Perkins  
**Client PO:** 141391

## Instrument Calibration

**Sequence:** 1100569  
**Instrument:** THG-05  
**Date:** 08/22/2011  
**Analyte:** Hg

**Total Mercury and Mercury Speciation by CVAFS**  
**Method:** EPA 1631


Lab ID	True Value	Result	Units	REC & Limits	
1100569-IBL1		6.42	pg of Hg		
1100569-IBL2		8.76	pg of Hg		
1100569-IBL3		10.70	pg of Hg		
1100569-IBL4		11.08	pg of Hg		
1100569-CAL1	25.00	25.81	pg of Hg	103%	
1100569-CAL2	100.0	92.73	pg of Hg	93%	
1100569-CAL3	500.0	517.2	pg of Hg	103%	
1100569-CAL4	2500	2552	pg of Hg	102%	
1100569-CAL5	10000	9933	pg of Hg	99%	
1100569-ICV1	1568	1468	pg of Hg	94%	85-115
1100569-CCB1		15.8	pg of Hg		
1100569-CCV1	500.0	524.8	pg of Hg	105%	77-123
1100569-CCV2	500.0	514.0	pg of Hg	103%	77-123



## Sample Containers

<b>Lab ID:</b> 1134009-01			<b>Report Matrix:</b> FGD Wastewater			<b>Collected:</b> 08/10/2011	
<b>Sample:</b> BioReactor 1 Inf			<b>Sample Type:</b> Sample			<b>Received:</b> 08/16/2011	
<b>Des</b>	<b>Container</b>	<b>Size</b>	<b>Lot</b>	<b>Preservation</b>	<b>P-Lot</b>	<b>pH</b>	<b>Ship. Cont.</b>
A	Bottle FLPE Hg-T	125 mL	71313080 60	None	N/A		Cardboard Box
<b>Lab ID:</b> 1134009-02			<b>Report Matrix:</b> DIW			<b>Collected:</b> 08/02/2011	
<b>Sample:</b> Hg Blk BioReactor 1 Inf			<b>Sample Type:</b> Field Blank			<b>Received:</b> 08/16/2011	
<b>Des</b>	<b>Container</b>	<b>Size</b>	<b>Lot</b>	<b>Preservation</b>	<b>P-Lot</b>	<b>pH</b>	<b>Ship. Cont.</b>
A	Bottle FLPE Hg-T	125 mL	71313080 60	0.1% HCl (BRL)	1121032	<2	Cardboard Box
<b>Lab ID:</b> 1134009-03			<b>Report Matrix:</b> FGD Wastewater			<b>Collected:</b> 08/10/2011	
<b>Sample:</b> BioReactor 2 Inf			<b>Sample Type:</b> Sample			<b>Received:</b> 08/16/2011	
<b>Des</b>	<b>Container</b>	<b>Size</b>	<b>Lot</b>	<b>Preservation</b>	<b>P-Lot</b>	<b>pH</b>	<b>Ship. Cont.</b>
A	Bottle FLPE Hg-T	250 mL	71313080 60	None	N/A		Cardboard Box
<b>Lab ID:</b> 1134009-04			<b>Report Matrix:</b> DIW			<b>Collected:</b> 08/02/2011	
<b>Sample:</b> Hg Blk BioReactor 2 Inf			<b>Sample Type:</b> Field Blank			<b>Received:</b> 08/16/2011	
<b>Des</b>	<b>Container</b>	<b>Size</b>	<b>Lot</b>	<b>Preservation</b>	<b>P-Lot</b>	<b>pH</b>	<b>Ship. Cont.</b>
A	Bottle FLPE Hg-T	125 mL	71313080 60	0.1% HCl (BRL)	1121032	<2	Cardboard Box
<b>Lab ID:</b> 1134009-05			<b>Report Matrix:</b> FGD Wastewater			<b>Collected:</b> 08/10/2011	
<b>Sample:</b> BioReactor 2 Eff			<b>Sample Type:</b> QC Sample			<b>Received:</b> 08/16/2011	
<b>Des</b>	<b>Container</b>	<b>Size</b>	<b>Lot</b>	<b>Preservation</b>	<b>P-Lot</b>	<b>pH</b>	<b>Ship. Cont.</b>
A	Bottle FLPE Hg-T	250 mL	71313080 60	None	N/A		Cardboard Box
<b>Lab ID:</b> 1134009-06			<b>Report Matrix:</b> DIW			<b>Collected:</b> 08/02/2011	
<b>Sample:</b> Hg Blk BioReactor 2 Eff			<b>Sample Type:</b> Field Blank			<b>Received:</b> 08/16/2011	
<b>Des</b>	<b>Container</b>	<b>Size</b>	<b>Lot</b>	<b>Preservation</b>	<b>P-Lot</b>	<b>pH</b>	<b>Ship. Cont.</b>
A	Bottle FLPE Hg-T	125 mL	71313080 60	0.1% HCl (BRL)	1121032	<2	Cardboard Box

**Project ID:** DUK-HV1101  
**PM:** Tiffany Stilwater

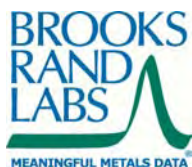

Analytical Lab  
Page 24 of 29  
**Client PM:** Jay Perkins  
**Client PO:** 141391

## Shipping Containers

### Cardboard Box

**Received:** August 16, 2011 9:00  
**Tracking No:** 4726 7966 2932 via FedEx  
**Coolant Type:** None  
**Temperature:** ambient

**Description:** Cardboard Box  
**Damaged in transit?** No  
**Returned to client?** No

**Custody seals present?** No  
**Custody seals intact?** No  
**COC present?** Yes



## CHAIN OF CUSTODY RECORD AND ANALYSIS REQUEST FORM

1134009



Duke Energy Analytical Laboratory

Mail Code MGO342 (Building 7405)

13339 Hagers Ferry Rd  
Huntersville, N.C. 28078

(704) 875-4245

Fax: (704) 875-4349

Analytical Lab  
Page 25 of 29

Customer must complete

1) Project Name	Bellevue - FGD	2) Phone No:
3) Client:	WWTS (2011, Bi-Weekly Sampling) Bill Kennedy, Melonie Martin, Wayne Chapman, Tom Johnson *	4) Fax No:
5) Business Unit:		Mail Code:
6) Oper. Unit:	9) Res. Type:	10) Reso. Center:

LAB USE ONLY
1) Lab ID

Customer to complete appropriate columns to right

Se Specification Bottle	ID	13 Sample Description or ID	Date	Time	Signature	17 Comp.	18 Grab	Hg 1631 (sample 2nd week only)
		BioReactor 1 Inf	8-10-11	1220	[Signature]			1
		Hg Blk BioReactor 1 Inf	8/10/11	1415	[Signature]			1
		BioReactor 2 Inf	8/10/11	1230	[Signature]			1
		Hg Blk BioReactor 2 Inf	8/10/11	1415	[Signature]			1
		BioReactor 2 Eff	8-10-11	1225	[Signature]			1
		Hg Blk BioReactor 2 Eff	8/10/11	1415	[Signature]			1

Use the BioReactor 2 Inf or EFF sample as the MS/MSD

1) Relinquished By	8-10-11	1525	2) Accepted By	8-10-11	1525	22) Requested Turnaround
3) Relinquished By			4) Accepted By	8/10/11	0900	14 Days
5) Relinquished By	15	0900	6) Accepted By			7) Days
7) Relinquished By	8-11-11	1300	8) Accepted By			48 Hr
9) Seal/Locked By	8-11-11		10) Seal/Lock Opened By			*Other
11) Seal/Locked By	15		12) Seal/Lock Opened By			*Add. Cost Will Apply

Comments: \* Metals=As, Ag, B, Cu, Cr, Ni, Se, Zn \*thomas.d.johnson@siemens.com

Customer, IMPORTANT!  
Please indicate desired turnaround.

Analytical Laboratory Use Only

Sample Class OTHER

ORDER # 511080089

Logged By [Signature]

Date &amp; Time 8-10-11 1540

Val

Brooks Rand

PO PO141391

Cooler Temp (C)

Preserv.: 1=HCL, 2=H2SO4, 3=HNO3, 4=Ice, 5=None

MR #

Customer to complete all appropriate non-shaded areas.

16 Analyses Required

Samples - Contained 2nd Wednesday each month

Date

Time

Signature

17 Comp.

18 Grab

HCL

H2SO4

HNO3

Ice

None

NC

SC

SAMPLE PROGRAM

Water

Drinking Water

RORA Waste

UST

Ground

NPDES

5

ORIGINAL TO LAB

COPY TO CLIENT

Page 2 of 2

DISTRIBUTION

ORIGINAL TO LAB

COPY TO CLIENT

140-334 N-1 01-08 00

ND BFI A

98107

MA-US

SEA

4726 7966 2932

TRK# 0201

PRIORITY OVERNIGHT

TUE - 16 AUG A1

020105050101024

E

FedEx

Express

020105050101024

020105050101024

SHIP DATE: 15AUG11

ACTWGT: 3.1 LB

CAD: 798987/CRFE2472

DIMS: 13x9x6 IN

BILL SENDER

ORIGIN ID: SRWA (980) 875-5213

DUKE ENERGY

G.C. SHARMA

13339 HAGERFERRY RD

BLDG # 7405

HUNTERSVILLE, NC 28078

UNITED STATES US

ATTN: MICHELLE BRISCOE

BROOKS RAND

3958 6TH AVENUE NW

SEATTLE WA 98107

REF: (206) 632-6206

DEPT: INV: 1

**From:** Perkins, Jay C  
**To:** [Chapman, Wayne C](#)  
**Subject:** Br  
**Date:** Monday, August 15, 2011 1:55:00 PM  
**Attachments:** [image001.gif](#)

---

FGD Purge EFF 97 mg/L  
BioReactor 2 Eff 95 mg/L

Sampled: August 10

Wayne, this is preliminary data.

Thanks,

Jay Perkins  
Duke Energy Analytical Laboratory  
W-980-875-5348  
C-704-796-6598



# CHAIN OF CUSTODY RECORD AND ANALYSIS REQUEST FORM

Analytical Lab  
Page 28 of 29



**Duke Energy Analytical Laboratory**  
Mail Code MGO3A2 (Building 7405)  
13339 Hagers Ferry Rd  
Huntersville, N. C. 28078  
(704) 875-5245  
Fax: (704) 875-4349

## Analytical Laboratory Use Only

ORDER# <b>511080089</b>	MATRIX: OTHER	Samples Originating From	NC SC
Logged By <b>Cpk</b>	Date & Time <b>8-10-11 1546</b>	SAMPLE PROGRAM Water _____ Ground NPDES _____ Drinking Water UST _____ RCRA Waste _____	
Cooler Temp (C) <b>&lt; 1</b>			

<sup>19</sup>Page 1 of 2  
**DISTRIBUTION**  
ORIGINAL to LAB,  
COPY to CLIENT

1)Project Name <b>Belews - FGD WWTS Bi-Monthly Sampling)</b>	2)Phone No:
2) Client: <b>Bill Kennedy, Melonie Martin, Wayne Chapman, Tom Johnson **</b>	4)Fax No:
5)Business Unit:	6)Process: Mail Code:
8)Oper. Unit:	10)Reso. Center:

MR #		Customer to complete all appropriate non-shaded areas.		15 Analyses Required		17 Comp.		18 Grab		TDS		Hg - 245.1		Br (Dionex)		Metals*		Se, soluble (no dig.)		Se, speciation - vendor to AS&C (Important to place filled bottle back into both baggies)	
		Sampling conducted: 2nd and 4th Wednesday																			
Date	Time	Signature																			
8/10	0800	Tom Johnson																			
8/10	0800	Tom Johnson																			
8/10	0800																				
8/10	0800																				
8/10	0800																				
8/2/11	1415	R. Davis																			
8/2/11	1415	R. Davis																			

Filtering of the Se is performed in the field please provide a filter blank too.

## LAB USE ONLY

<sup>11</sup>Lab ID

2011016915  
16  
17  
18  
19  
20  
21

Customer to complete appropriate columns to right

Se Speciation Bottle ID	<sup>13</sup> Sample Description or ID	Date	Time	Signature
B10871	FGD Purge Eff	8/10	0800	Tom Johnson
	EQ Tank Eff.	8/10	0800	Tom Johnson
B10859	BioReactor 1 Inf	8/10	0800	
	BioReactor 2 Inf	8/10	0800	
B11239	BioReactor 2 Eff	8/10	0800	
	Filter Blk	8/2/11	1415	R. Davis
B11716	Metals Trip Blk	8/2/11	1415	R. Davis

Customer to sign & date below - fill out from left to right.

1) Relinquished By <i>Tom Johnson</i>	Date/Time <i>8/10/11 0900</i>	2) Accepted By <i>Tom Johnson</i>	Date/Time <i>8-10-11 1215</i>
3) Relinquished By <i>Tom Johnson</i>	Date/Time <i>8-10-11 1525</i>	4) Accepted By <i>Candy Knox</i>	Date/Time <i>8-10-11 1525</i>
5) Relinquished By <i>Cpk</i>	Date/Time <i>8-11-11 1300</i>	6) Accepted By <i>Cpk</i>	Date/Time <i>8-11-11 1300</i>
7) Relinquished By <i>Cpk</i>	Date/Time <i>8-11-11 1300</i>	8) Accepted By <i>Cpk</i>	Date/Time <i>8-11-11 1300</i>
9) Seal/Locked By <i>Cpk</i>	Date/Time <i>8-11-11 1300</i>	10) Seal/Lock Opened By <i>Cpk</i>	Date/Time <i>8-11-11 1300</i>
11) Seal/Locked By <i>Cpk</i>	Date/Time <i>8-11-11 1300</i>	12) Seal/Lock Opened By <i>Cpk</i>	Date/Time <i>8-11-11 1300</i>

Comments

\* B by ICP As, Cr, Cu, Ni, Se, Ag, Zn by IMS Digestions = TRM thomas.d.johnson@siemens.com

Customer, IMPORTANT!  
Please indicate desired turnaround.

## <sup>22</sup>Requested Turnaround

14 Days \_\_\_\_\_

\*7 Days \_\_\_\_\_

\* 48 Hr \_\_\_\_\_

\*Other \_\_\_\_\_

\* Add. Cost Will Apply



# CHAIN OF CUSTODY RECORD AND ANALYSIS REQUEST FORM

Analytical Lab  
Page 29 of 29



**Duke Energy Analytical Laboratory**  
Mail Code MGO3A2 (Building 7405)  
13339 Hagers Ferry Rd  
Huntersville, N. C. 28078  
(704) 875-5245  
Fax: (704) 875-4349

## Analytical Laboratory Use Only

ORDER # <b>J11080089</b>	Sample Class <b>OTHER</b>	Samples Originating From NC _____ SC _____
Logged By <b>cpt</b>	Date & Time <b>8-10-11 1546</b>	SAMPLE PROGRAM Water _____ Drinking Water UST _____ RCRA Waste _____
Ver <b>Brooks Rand</b>	PO <b>PO#141391</b>	Cooler Temp (C) <b>&lt;1</b> Preserv.: 1=HCL 2=H <sub>2</sub> SO <sub>4</sub> 3=HNO <sub>3</sub> 4=Ice 5=None

<sup>19</sup>Page 2 of 2  
**DISTRIBUTION**  
ORIGINAL to LAB,  
COPY to CLIENT

1) Project Name <b>Belews - FGD</b> <b>WWTS (2011, Bi-Weekly Sampling)</b>	2) Phone No:
2) Client: <b>Bill Kennedy, Melonie Martin, Wayne Chapman, Tom Johnson *</b>	4) Fax No:
5) Business Unit:	6) Process: Mail Code:
8) Oper. Unit:	9) Res. Type: 10) Reso. Center:

MR #						16 Analyses Required													
Customer to complete all appropriate non-shaded areas.																			
Sampling conducted: 2nd Wednesday each month																			
Se Speciation Bottle ID		13 Sample Description or ID		Date	Time	Signature	17 Comp.	18 Grab											
		BioReactor 1 Inf		8-10-11	1220	D. Morris													
		Hg Blk BioReactor 1 Inf		8/2/11	1415	R. Davis													
		BioReactor 2 Inf		8-10-11	1230	D. Morris													
		Hg Blk BioReactor 2 Inf		8/2/11	1415	R. Davis													
		BioReactor 2 Eff		8-10-11	1225	D. Morris													
		Hg Blk BioReactor 2 Eff		8/2/11	1415	R. Davis													
Use the Bioreactor 2 Inf or EFF sample as the MS/MSD																			

LAB USE ONLY  
11 Lab ID  
**2011016934**  
**35**  
**36**  
**37**  
**38**  
**39**

Customer to complete appropriate columns to right

Customer to sign & date below - fill out from left to right.

1) Relinquished By <b>D. Morris</b>	Date/Time <b>8-10-11 1525</b>	2) Accepted By <b>Cindy Knox</b>	Date/Time <b>8-10-11 1525</b>
3) Relinquished By	Date/Time	4) Accepted By	Date/Time
5) Relinquished By <b>15 cpt</b>	Date/Time	6) Accepted By:	Date/Time
7) Relinquished By <b>cpt/ld</b>	Date/Time <b>8-11-11 1300</b>	8) Accepted By:	Date/Time
9) Seal/Locked By <b>cpt/ld</b>	Date/Time <b>8-11-11</b>	10) Seal/Lock Opened By	Date/Time
11) Seal/Locked By <b>15</b>	Date/Time	12) Seal/Lock Opened By	Date/Time

**Customer, IMPORTANT!**  
Please indicate desired turnaround.

**22 Requested Turnaround**

14 Days \_\_\_\_\_

\*7 Days  
**8-23-11**

\*48 Hr \_\_\_\_\_

\*Other \_\_\_\_\_  
\* Add. Cost Will Apply

Comments  
\* Metals=As, Ag, B, Cu, Cr, Ni, Se, Zn \*thomas.d.johnson@siemens.com