## Supporting information

Legacy and Emerging Perfluoroalkyl Substances Are Important

Drinking Water Contaminants in the Cape Fear River Watershed

of North Carolina

Supporting information includes analytical method description, 6 tables, and 5 figures.

Mei Sun<sup>1,2,\*</sup>, Elisa Arevalo<sup>2</sup>, Mark Strynar<sup>3</sup>, Andrew Lindstrom<sup>3</sup>, Michael Richardson<sup>4</sup>, Ben Kearns<sup>4</sup>, Adam Pickett<sup>5</sup>, Chris Smith<sup>6</sup>, and Detlef R.U. Knappe<sup>2</sup>

- <sup>1</sup>Department of Civil and Environmental Engineering University of North Carolina at Charlotte Charlotte, North Carolina 28223, USA
- <sup>2</sup> Department of Civil, Construction, and Environmental Engineering North Carolina State University Raleigh, North Carolina 27695, USA
- National Exposure Research Laboratory
   U.S. Environmental Protection Agency
   Research Triangle Park, North Carolina 27711, USA
- <sup>4</sup> Cape Fear Public Utility Authority Wilmington, North Carolina 28403, USA
- <sup>5</sup> Town of Pittsboro Pittsboro, North Carolina 27312, USA
- <sup>6</sup> Fayetteville Public Works Commission Fayetteville, North Carolina 28301, USA

<sup>\*</sup>Corresponding Author Email: msun8@uncc.edu; Phone: 704-687-1723

Analytical standards: PFASs studied in this research are listed in Table S1. For legacy PFASs, native and isotopically labeled standards were purchased from Wellington Laboratories (Guelph, Ontario, Canada). Native PFPrOPrA was purchased from Thermo Fisher Scientific (Waltham, MA). No analytical standards were available for other PFECAs.

PFAS quantification: PFAS concentrations in samples from DWTPs and adsorption tests were determined by liquid chromatography tandem mass spectrometry (LC-MS/MS) using a large-volume (0.9 mL) direct injection method. An Agilent 1100 Series LC pump and PE Sciex API 3000 LC-MS/MS system equipped with a 4.6 mm x 50 mm HPLC column (Kinetex C18 5μm 100Å, Phenomenex Inc.) was used for PFAS analysis. The eluent gradient is shown in Table S4 in SI. All samples, calibration standards, and quality control samples were spiked with isotopically labeled internal standards, filtered through 0.45-μm glass microfiber syringe filters, and analyzed in duplicate. The MS transitions for PFAS analytes and internal standards are shown in Table S5 in SI. The quantitation limit (QL) was 25 ng/L for PFOS and perfluorodecanoic acid, and 10 ng/L for other legacy PFASs and PFPrOPrA. The QL was defined as the first point of the standard curve, for which the regression equation yielded a calculated value within ±30% error. For PFECAs without analytical standards, chromatographic peak areas are reported.

PFAS concentrations along the treatment train of DWTP C were analyzed using a Waters Acquity ultra performance liquid chromatograph interfaced with a Waters Quattro Premier XE triple quadrupole mass spectrometer (Waters, Milford, MA, USA) after solid phase extraction. Method details are described elsewhere. The QL for all PFASs with analytical standards was 0.2 ng/L, and peak areas were recorded for PFECAs without standards.

Table S1. Perfluoroalkyl substances (PFASs) detected in the Cape Fear River (CFR) watershed

Compound	Molecular weight	Formula	CAS#	# of perfluorinated carbons	Chain length (including all C, O and S)					
Perfluorocarboxylic acids (PFCAs)										
Perfluorobutanoic acid (PFBA)	214.0	C4HF7O2	375-22-4	3	4					
Perfluoropentanoic acid (PFPeA)	264.0	C5HF9O2	2706-90-3	4	5					
Perfluorohexanoic acid (PFHxA)	314.1	C <sub>6</sub> HF <sub>11</sub> O <sub>2</sub>	307-24-4	5	6					
Perfluoroheptanoic acid (PFHpA)	364.1	C7HF13O2	375-85-9	6	7					
Perfluorooctanoic acid (PFOA)	414.1	C <sub>8</sub> HF <sub>15</sub> O <sub>2</sub>	335-67-1	7	8					
Perfluorononanoic acid (PFNA)	464.1	C9HF17O2	375-95-1	8	9					
Perfluorodecanoic acid (PFDA)	514.1	C10HF19O2	335-76-2	9	10					
Perflu	orosulfonic a	cids (PFSAs)								
Perfluorobutane sulfonic acid (PFBS)	300.1	C <sub>4</sub> HF <sub>9</sub> SO <sub>3</sub>	375-73-5	4	5					
Perfluorohexane sulfonic acid (PFHxS)	400.1	C <sub>6</sub> HF <sub>13</sub> SO <sub>3</sub>	355-46-4	6	7					
Perfluorooctane sulfonic acid (PFOS)	500.1	C8HF17SO3	1763-23-1	8	9					
Perfluoroalkyl ether carboxyli	c acids with o	ne ether grou	p (mono-ether	PFECAs)						
Perfluoro-2-methoxyacetic acid (PFMOAA)	180.0	C3HF5O3	674-13-5	2	4					
Perfluoro-3-methoxypropanoic acid (PFMOPrA)	230.0	C <sub>4</sub> HF <sub>7</sub> O <sub>3</sub>	377-73-1	3	5					
Perfluoro-4-methoxybutanoic acid (PFMOBA)	280.0	C5HF9O3	863090-89-5	4	6					
Perfluoro-2-propoxypropanoic acid (PFPrOPrA)	330.1	C6HF11O3	13252-13-6	5	7					
Perfluoroalkyl ether carboxylic acids with multiple ether group (multi-ether PFECAs)										
Perfluoro(3,5-dioxahexanoic) acid (PFO2HxA)	246.0	C4HF7O4	39492-88-1	3	6					
Perfluoro(3,5,7-trioxaoctanoic) acid (PFO3OA)	312.0	C5HF9O5	39492-89-2	4	8					
Perfluoro(3,5,7,9-tetraoxadecanoic) acid (PFO4DA)	378.1	C <sub>6</sub> HF <sub>11</sub> O <sub>6</sub>	39492-90-5	5	10					

Table S2. Operational conditions of DWTP C on sampling day (August 18, 2014)

Parameter	Value
Raw water ozone dose	3.1 mg/L
Raw water total organic carbon concentration	6.0 mg/L
Aluminum sulfate coagulant dose	43 mg/L
Coagulation pH	5.70
Settled water ozone dose	1.3 mg/L
Settled water total organic carbon concentration	1.90 mg/L
Empty bed contact time in biological activated carbon filters	9.4 minutes for granular activated carbon layer 2.3 minutes for sand layer
Medium pressure UV dose	25 mJ/cm <sup>2</sup>
Free chlorine dose	1.26 mg/L as Cl <sub>2</sub>
Free chlorine contact time	17.2 hours

Table S3. Water quality characteristics of surface water used in adsorption tests

Non-purgeable	Ultraviolet absorbance		Alkalinity	Conductivity	
organic carbon (mg/L)	at a wavelength of 254 nm		(mg/L as CaCO <sub>3</sub> )	(µS/cm)	
9.036	0.399	7.53	19	133.5	

Table S4. LC gradient method for PFAS analysis

Time (min)	Mobile Phase A% (v/v)	Mobile Phase B%	Flow Rate (mL/min)
0 – 2	95	5	0.9
2 – 5	95	5	0.9
5 – 10	95 → 10	$5 \rightarrow 90$	0.9
10 – 10.1	10	90	0.9
10.1 – 14	$10 \rightarrow 95$	$90 \rightarrow 5$	0.9

Mobile phase A: 2 mM ammonium acetate in ultrapure water with 5% methanol

Mobile phase B: 2 mM ammonium acetate in acetonitrile with 5% ultrapure water

Table S5. MS transitions for PFAS Analysis

	Compound	MS/MS Transition	Internal standard	
	PFBA	$212.8 \rightarrow 168.8$	13C4-PFBA	
	PFPeA	$262.9 \rightarrow 218.8$	13C2- PFHxA	
	PFHxA	$313.6 \rightarrow 268.8$	13C2- PFHxA	
Legacy PFASs PFECAs	РҒНрА	$362.9 \rightarrow 318.8$	13C4- PFOA	
L DEAC-	PFOA	$413.0 \rightarrow 368.8$	13C4- PFOA	
Legacy PFASs	PFNA	$463.0 \to 418.8$	13C4- PFOA	
	PFDA	$513.1 \rightarrow 68.8$	13C2-PFDA	
PFBS PFHxS	PFBS	$299.1 \to 98.8$	18O2-PFHxS	
	PFHxS	$399.1 \to 98.8$	18O2-PFHxS	
	PFOS	$498.9 \to 98.8$	13C4-PFOS	
	PFMOAA	$180.0 \to 85.0$	N/A	
	PFMOPrA	$229.1 \rightarrow 184.9$	N/A	
	PFMOBA	$279.0 \rightarrow 234.8$	N/A	
PFECAs	PFPrOPrA	$329.0 \rightarrow 284.7$	13C2- PFHxA	
	PFO2HxA	$245.1 \to 85.0$	N/A	
	PFO3OA	311. → 84.9	N/A	
	PFO4DA	$377.1 \to 85.0$	N/A	
	Perfluoro-n-[1,2,3,4-13C4]butanoic acid (13C4-PFBA)			
Internal standards	Perfluoro-n-[1,2- <sup>13</sup> C <sub>2</sub> ]hexanoic acid (13C2-PFHxA)	$315.1 \rightarrow 269.8$		
	Perfluoro-n-[1,2,3,4- <sup>13</sup> C <sub>2</sub> ]octanoic acid (13C4-PFOA)	$417.0 \rightarrow 372.0$	NI-t - multi-ski	
	Perfluoro-n-[1,2- <sup>13</sup> C <sub>2</sub> ]decanoic acid (13C2-PFDA)	515.1 → 469.8	Not applicable	
	Sodium perfluoro-1- hexane[18O2]sulfonate (18O2-PFHxS)	403.1 → 83.8		
	Sodium perfluoro-1-[1,2,3,4- <sup>13</sup> C <sub>4</sub> ]octane sulfonate (13C4-PFOS)	502.9 → 79.9		

Table S6. Maximum, minimum, mean and median concentrations (ng/L) of PFASs at three drinking water intakes. \*

		Com	munity A		Community B			Community C				
	max	min	median	mean	max	min	median	mean	max	min	median	mean
PFBA	99	<10	26	33	38	<10	12	12	104	<10	12	22
PFPeA	191	14	44	62	38	<10	19	19	116	<10	30	36
PFHxA	318	<10	48	78	42	<10	<10	11	24	<10	<10	<10
PFHpA	324	<10	39	67	85	<10	<10	11	24	<10	<10	<10
PFOA	137	<10	34	46	32	<10	<10	<10	17	<10	<10	<10
PFNA	38	<10	<10	<10	<10	<10	<10	<10	<10	<10	<10	<10
PFDA	35	<25	<25	<25	<25	<25	<25	<25	<25	<25	<25	<25
PFBS	80	<10	<10	<10	11	<10	<10	<10	<10	<10	<10	<10
PFHxS	193	<10	10	14	14	<10	<10	<10	14	<10	<10	<10
PFOS	346	<25	29	44	43	<25	<25	<25	40	<25	<25	<25
PFPrOPrA	<10	<10	<10	<10	10	<10	<10	<10	4560	55	304	631
PFOA+PFOS	447	0	64	90	59	0	0	9	55	<10	<10	<10
∑ PFASs**	1502	18	212	355	189	0	47	62	4696	55	345	710

<sup>\*</sup> Concentrations less than quantitation limits were considered as zero to calculate means and  $\Sigma$  PFASs.

<sup>\*\*</sup> Other PFECAs were present in water samples from community C but could not be quantified and were therefore not included in  $\Sigma$  PFASs

Figure S1. Molecular structures of PFECAs evaluated in this study

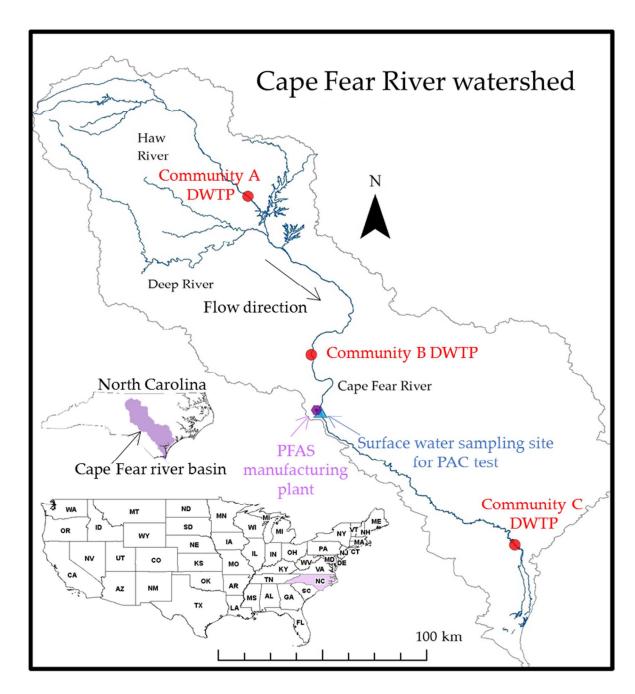
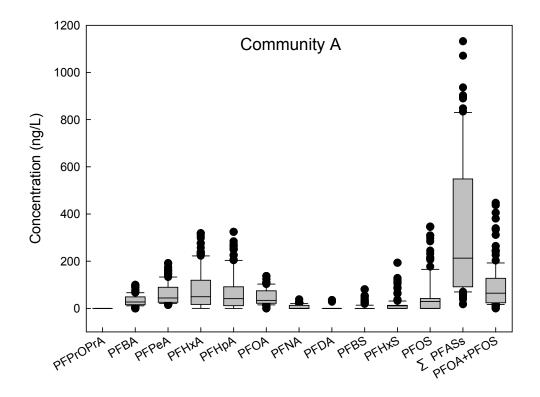
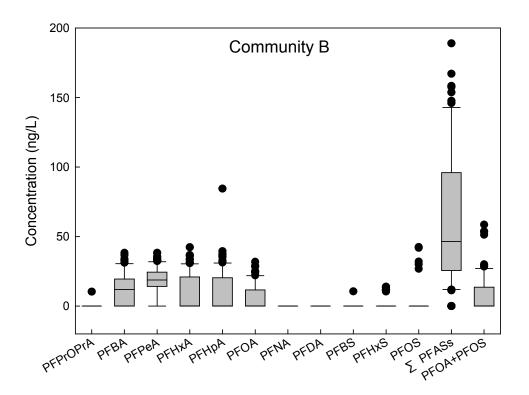


Figure S2. Sampling sites in the Cape Fear River watershed, North Carolina. The scale is for the Cape Fear River watershed map.





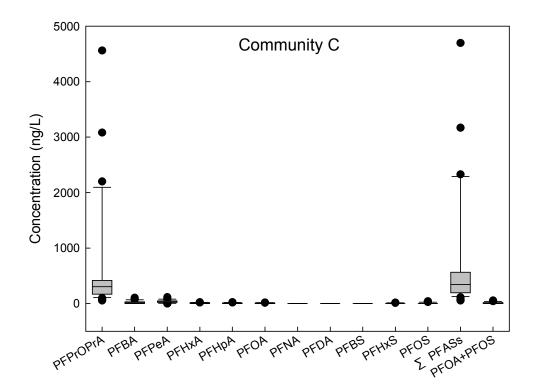


Figure S3. PFAS concentration distributions in the CFR watershed at three drinking water intakes. Concentrations less than quantitation limits were considered as zero. Upper and lower edges of a box represent the 75<sup>th</sup> and 25<sup>th</sup> percentile, respectively; the middle line represents the median; upper and lower bars represent the 90<sup>th</sup> and 10<sup>th</sup> percentile, respectively; and dots represent outliers (>90<sup>th</sup> or <10<sup>th</sup> percentile).

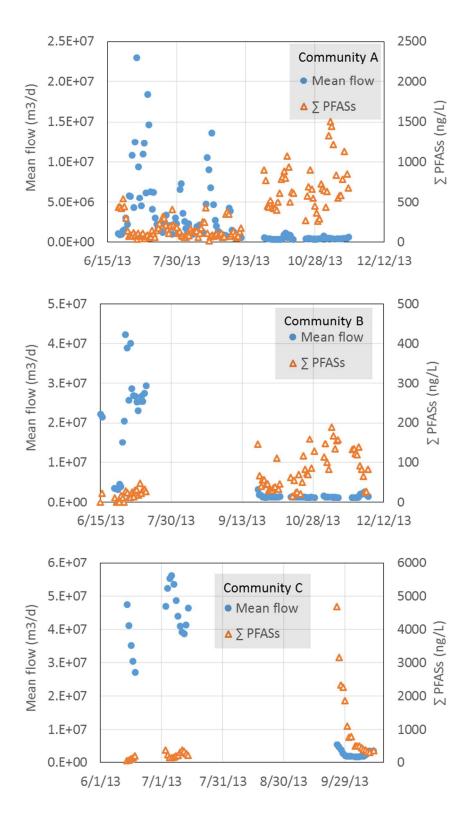


Figure S4. Total PFAS concentrations in the source water and stream flow at the three studied DWTPs. Stream flow data were acquired from US Geological Survey stream gage records

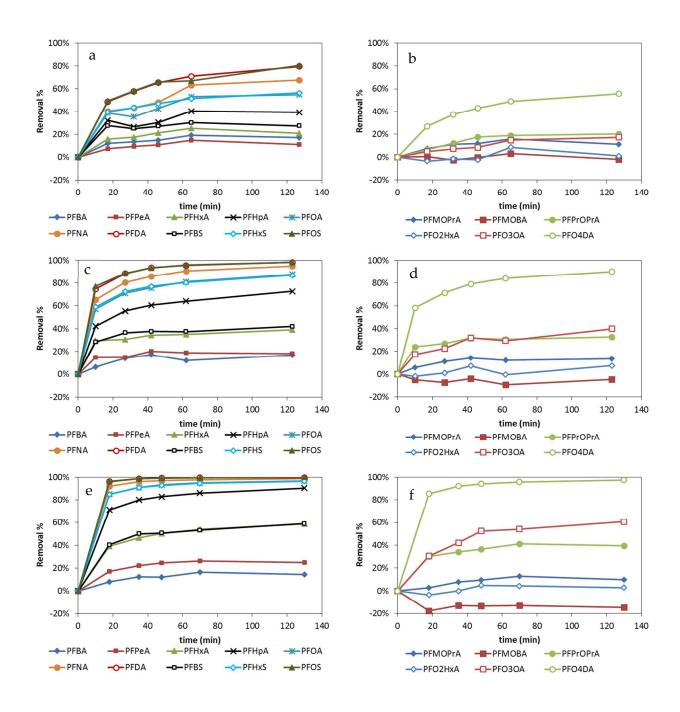


Figure S5. PFAS adsorption at powdered activated carbon doses of (a, b) 30 mg/L, (c, d) 60 mg/L and (e, f) 100 mg/L. Figures show average PFAS removal percentages of duplicate tests.

## Reference

1. Nakayama, S.; Strynar, M. J.; Helfant, L.; Egeghy, P.; Ye, X.; Lindstrom, A. B., Perfluorinated compounds in the Cape Fear drainage basin in North Carolina. *Environ. Sci. Technol.* **2007**, *41*, (15), 5271-5276.