#### UNITED STATES ENVIRONMENTAL PROTECTION AGENCY WASHINGTON, D.C. 20460



OFFICE OF CHEMICAL SAFETY AND POLLUTION PREVENTION

#### **MEMORANDUM**

Date: 6/29/2015

SUBJECT: 2,4-D: Data Evaluation Records (DERs) for EDSP Tier 1 Assays

PC Code: 030001 DP Barcode: D398637, D398638, D398640

Decision No.: NA

Petition No.: NA

Registration No.: NA

Regulatory Action: NA

Risk Assessment Type: NA
TXR No.: 0052104
CAS No.: 94-75-7
MRID No.: See Table
Case No.: NA
CAS No.: 94-75-7
40 CFR: NA

Ver.Apr. 2010

FROM: Greg Akerman, Ph.D.

Immediate Office

Health Effects Division (7509P)

THROUGH: Jess Rowland Jess Row

Deputy Director

Health Effects Division

TO: Jolene Trujillo

Biologist/Chemical Review Manager

Risk Management and Implementation Branch V

Pesticide Re-evaluation Division (7505P)

#### I. ACTION REQUESTED

The Pesticide Re-evaluation Division (PRD) of OPP has requested that the Health Effects Division (HED) review the Endocrine Disruptor Screening Program (EDSP) Tier 1 assays submitted in response to the agency's Test Order for 2,4-D: Test Order # CON-030001-1.

#### II. RESPONSE

Attached are the EDSP Tier 1 assay DERs for 2,4-D.

#### III. MRID Table

Chemical:	2,4-D	PC Code: 030001
Guideline	Assay	MRID
890.1100	Amphibian Metamorphosis Assay (Frog)	48317002
890.1150	Androgen Receptor Binding (Rat Prostate)	48614301
890.1200	Aromatase Assay (Human Recombinant)	48614302
890.1250	Estrogen Receptor Binding	48614303
890.1300	Estrogen Receptor Transcriptional Activation (Human Cell Line HeLa-9903)	48614304
890.1350	Fish Short-Term Reproduction	48317001
890.1400	Hershberger (Rat)	NA
890.1450	Female Pubertal (Rat)	NA
890.1500	Male Pubertal (Rat)	NA
890.1550	Steroidogenesis (Human Cell Line – H295R)	48614305
890.1600	Uterotrophic (Rat)	NA

NA= Not Applicable (Requirement satisfied by other scientifically relevant information (OSRI)).

Digitally signed by ROBIN STERNBERG DN: c=US, o=U.S. Government, ou=USEPA, ou=Staff, cn=ROBIN STERNBERG, dnQualifier=0000039126 Date: 2015.06.01 13:09:51-04'00'

Data Requirement: EPA DP Barcode 388579

OECD Data Point 231

EPA MRID 48317002 EPA Guideline 890.1100

Amphibian Metamorphosis Assay (Frog)

Test Material: 2,4-Dichlorophenoxyacetic acid Purity (%): 98.6

Common Name 2,4-D
Chemical Name IUPAC

CAS Name

CAS No. 94-75-7 Synonyms 2,4-D

2,4-D acid

EPA PC Code 030001

Primary Reviewer: Catherine Aubee

Signature:

2015.06.01

13:15:13

-04'00'

Biologist, USEPA/OCSPP/OPP/EFED/ERB1 Date: 07/30/2012

Additional Reviewer: Alicia Korol Signature: No longer with EPA

Biologist, USEPA/OCSPP/OPP/EFED/ERB1 Date: 06/24/2011

Additional Reviewer: Anita Ullagaddi Signature: No longer with EPA

Biologist, USEPA/OCSPP/OPP/EFED/ERB1 Date: 02/24/2012

Final Additional Reviewer: Robin Sternberg Signature:

Wildlife Biologist, USEPA/OCSPP/OPP/EFED/ERB1 Date: 05/27/2015

Date Evaluation Completed: 05/27/2015

Page 1 of 61

DER Template Version: 22 September 2011

CITATION: Coady, K.K., T.A. Marino, and J. Thomas. November 29, 2010. 2,4-Dichlorophenoxyacetic Acid: The Amphibian Metamorphosis Assay Using the South African Clawed Frog, Xenopus laevis. The Dow Chemical Company. Midland, Michigan. Laboratory Project Study ID 101025. Sponsor: Industry Task Force II on 2,4-D Research Data.

The US EPA Endocrine Disruptor Screening Program (EDSP) Tier 1 screening battery is comprised of eleven screening assays intended to identify a chemical's likely endocrine bioactivity, i.e., its potential to interact with the estrogen, androgen, or thyroid (E, A, or T) pathways. The robustness of the Tier 1 battery is based on the strengths of each individual assay to identify potential endocrine bioactivity with complementary endpoints within the assay, where available, and redundancy across the battery. Thus, the results of each individual assay should not be considered in isolation but rather should be considered in the context of other assays in the battery as well as Other Scientifically Relevant Information (OSRI). In order to determine if a chemical has the potential to interact with the E, A or T pathways, a Weight of Evidence (WoE) evaluation of Tier 1 assay results, in combination with the findings in the OSRI, should be undertaken (refer to the WoE Document).

Disclaimer: The guideline recommendations in this DER template are offered as a general reference to aid in preparation of the DER. The purpose of these recommendations is not to serve as substitute for the Test Guidelines, nor to provide any guidance on how the study should be conducted.

Data Evaluation Record on the Toxicity of 2,4-D to Amphibians, Metamorphosis Assay

EPA MRID Number 48317002

**EXECUTIVE SUMMARY** 

The 21-day assay of 2,4-D (purity 98.6%) on amphibian metamorphosis of South African Clawed Frog

(Xenopus laevis) was conducted under flow-through conditions. Amphibian larvae at Nieuwkoop-Faber (NF)

stage 51 (80/control and treatment group; 20/replicate) were exposed to a negative control and test chemical

nominal concentrations of 0.4, 4.0, 40.0, and 100.0 mg a.i./L. Mean measured concentrations were <0.120

(<LOQ, negative control), 0.273, 3.24, 38.0, and 113 mg a.i./L. The test system was maintained at 21.9 to

22.7 °C and a pH of 7.0 to 7.8.

Only one incidence of tadpole mortality occurred in the mid-low treatment group; the cause of death was

unknown. No clinical signs of toxicity were noted.

2,4-D did not affect Day 7 normalized (for snout-vent length) hind-limb length (HLL). However, there was a

statistically significant (p<0.05) decrease of 15% in Day 21 normalized HLL at the highest treatment level

compared with the negative control. There was no significant effect on median NF developmental stage, snout-

vent length (SVL), or body weight at Day 7 or Day 21. Asynchronous development was not observed. There

were no effects on thyroid gland histopathology. Late stage (>NF stage 60) tadpoles were observed in the

negative control and in all treatment levels; consistent with the guideline recommendations, these tadpoles

were excluded from analyses of growth and normalized HLL.

The study met all performance and validity criteria with the exception that the coefficients of variation (CVs)

for measured concentrations of the low and low-mid treatment groups were 56 and 22%, respectively,

exceeding the guideline performance criterion of <20%. This was likely due to biodegradation of the test

material in the test vessels.

The assay satisfies the EDSP Tier 1 Test Order requirements for an Amphibian Metamorphosis Assay (OCSPP

Guideline 890.1100).

Page 3 of 61

DER Template Version: 22 September 2011

Page 5 of 221

#### Results Synopsis:

Test organism NF stage at test initiation: 51

Test organism total length at test initiation (optional): Not reported

Test type: Flow-through

Table 1: Summary of Developmental and Thyroid Pathology/Histopathology Effects<sup>1,2</sup> in the Amphibian Metamorphosis Assay (AMA) with 2,4-D.

Treatment (mg a.i./L)	NF Develo	•	Hind Len	Limb gth <sup>3</sup>		nronous opment	Thyroid Gross and Histopathology
[mean-measured]	Day 7	Day 21	Day 7	Day 21	Day 7	Day 21	Day 21
0.273	No	No	No	No	No	No	No
3.24	No	No	No	No	No	No	No
38.0	No	No	No	No	No	No	No
113	No	No	No	Yes	No	No	No

A "yes" indicates a significant difference based on comparison to the negative (clean water) control, unless otherwise specified.

The criteria for significance are described in the Reviewer's Analysis and Statistical Verification sections of the DER. Conclusions regarding histopathology may be heavily weighted by the expert opinion of a board-certified pathologist.

<sup>&</sup>lt;sup>3</sup> Hind-limb length is normalized to snout-vent length (SVL).

#### I. MATERIALS AND METHODS

Guideline Followed:

This study was conducted following guidelines outlined in United States Environmental Protection Agency (USEPA) 2009, Endocrine Disrupter Screening Program Test Guidelines OPPTS 890.1100: Amphibian Metamorphosis (Frog). EPA 740-C-09-002, October 2009 and Organization for Economic Cooperation and Development (OECD). 2009. OECD Guideline for the Testing of Chemicals: The Amphibian Metamorphosis Assay. OECD231 Adopted 7 September 2009. The following deviations was noted:

 The CVs for measured concentrations of the low and low-mid treatment groups were 56 and 22%, respectively, exceeding the guideline performance criterion of ≤20%. This was likely due to biodegradation of the test material in the test vessels.

This deviation did not impact the interpretation of the study.

Compliance:

Signed and dated GLP, Quality Assurance and Data Confidentiality statements were provided. This study was conducted in compliance with the following: European Community (EC) — European Parliament and Council Directive 2004/1 O/EC (O.J. No. L 50/44, 20/02/2004); Organisation for Economic Co-Operation and Development (OECD) — OECD Series on Principles of Good Laboratory Practice and Compliance Monitoring, Number 1. OECD Principles on Good Laboratory Practice (as revised in 1 997) ENV/MC/CHEM(98)17; and US Environmental Protection Agency -- FIFRA GLPs Title 40 CFR, Part 160 - Federal Insecticide, Fungicide and Rodenticide Act (FIFRA); Good Laboratory Practice Standards, Final Rule.

A. Test Material

2,4-Dichlorophenoxyacetic acid

Page 5 of 61

DER Template Version: 22 September 2011

CAS No 94-75-7

Description:

OECD recommends describing water solubility, melting/boiling point stability in water and light, pKa, Pow or Kow, vapor pressure of test compound, expiration date.

Lot No./Batch No.: 2006 2433 8006-USA

Purity: 98.6%

Impurities: None reported

Stability of Compound: The CVs for mean-measured concentrations were 55, 22, 16, and 17% for

the low, low-mid, high-mid, and high concentrations, respectively. This was likely due to biodegradation of the test material in the test vessels. Overall recoveries were 68, 81, 95, and 113% of nominal for the low, low-mid, high-

mid, and high concentrations, respectively.

Storage Conditions of

Test Chemicals: Not reported. Aqueous stock solutions held in amber glassware or covered

during storage to prevent photodegradation.

#### B. Test Organism

Table 2: General Information About the Test Species and Parental Care.

Parameter	Value(s)	Details or Remarks	Guideline Recommendations
Species common name:	South African Clawed		EPA recommends African clawed frog
	Frog		(Xenopus laevis). Western [Africa] clawed
Species scientific name:	Xenopus laevis		frog Silurana (Xenopus) tropicalis may be
			used as an alternate species'; however, a
Species strain (if stated):			list of all of the necessary protocol
			deviations to accommodate this species is
			recommended for inclusion in the study
			report. The guideline recommends that the
			performance criteria used to support the
			reliability of the test be identified.

Disruptor Screening Program Tier 1 Assays (OCSPP Test Guideline Series 890). March 3, 2011. Office of Chemical Safety and Pollution Prevention U.S. Environmental Protection Agency (EPA). (2011). Corrections and Clarifications on Technical Aspects of the Test Guidelines for the Endocrine (OCSPP), Washington, D.C. (http://www.epa.gov/endo/pubs/assayvalidation/clarificationdoc.pdf).

Page 7 of 61

DER Template Version: 22 September 2011

Data Evaluation Record on the Toxicity of 2,4-D to Amphibians, Metamorphosis Assay

Parameter	Value(s)	Details or Remarks	Guideline Recommendations
Were parents maintained as in-	Yes	Breeding pairs were ordered from	EPA recommends that larvae used in the
house stock?		Xenopus Express (Brooksville, Florida).	assay be derived from in-house adults.
		Adult male and female frogs were	
		injected with human chorionic	
		gonadotropin at least twelve hours	
		before the desired breeding events.	
Were parental acclimation	No	Not reported	
conditions same as definitive			
test?			
Acclimation period for parental	days		
frogs (if applicable):			
Details on parental feeding:	Not reported		
Details on parental health:	Not reported		

Page 8 of 61

DER Template Version: 22 September 2011

# Data Evaluation Record on the Toxicity of 2,4-D to Amphibians, Metamorphosis Assay

Table 3: Larval Selection and Care.

Parameter	Value(s)	Parameter Value(s) Details or Remarks	Guideline Recommendations
Best single spawn?	oN	Not reported	EPA and OECD recommend that the best 2
Number of spawns evaluated (if applicable):		Not reported	- 3 individual spawns, with a minimum of 1500 larvae/spawn, be evaluated to identify
Number of eggs sampled per spawn:		Not reported	ine best single spawn, and that the larvae selected for testing originate from the best single spawn (i.e., the spawns are not comixed)
NF stage at test initiation	51		EPA recommends that the definitive study
Age at test initiation:	18 days post- fertilization (dpf)	Test initiation began when tadpoles were 16 days old (discrepancy)	be initiated with larvae at Nieuwkoop - Faber (NF) developmental stage 51 (17 days post-fertilization).
Mean total length at test initiation (if reported):	mm	Not reported	
Range of total length at test initiation (if reported):	mm : mm	Not reported	

Page 9 of 61

DER Template Version: 22 September 2011

Data Evaluation Record on the Toxicity of 2,4-D to Amphibians, Metamorphosis Assay

Parameter	Value(s)	Details or Remarks	Guideline Recommendations
Was the optional size selection method used?	No	Not reported	
Details on larval selection:	Not reported		
Loading rate (rearing density):		Not reported for larvae; 10 tadpoles/L	EPA recommends that rearing density (loading rate) not exceed approximately 10
			larvae/L culturing system for flow-through
			systems or 4 tadpoles/L in static-renewal
			exposure systems.
Type of food:	Sera Micron®		EPA recommends Sera Micron® throughout
Source of food:	Sera North America		pre-exposure (after NF stage 45/46) and
lodide concentration in diet (if	53.7 ug/g		during the entire 21-d definitive study. If another diet is used, the study report should
known):			provide analysis of iodide content and
			potential contaminants, and the diet should
			demonstrate equal performance to Sera
			Micron®.

Page 10 of 61

DER Template Version: 22 September 2011

Data Evaluation Record on the Toxicity of 2,4-D to Amphibians, Metamorphosis Assay

Parameter	Value(s)	Details or Remarks	Guideline Recommendations
Frequency of feeding:	2 times/day		EPA recommends that feeding occur at
			least twice per day.
Details on feeding regime:	The feeding regimen		It is recommended that food rations during
	was the same as		the pre-exposure period be increased along
	guideline		with larval growth to approximately 30
	recommendations		mg/larva/day by test initiation. EPA and
	(Table 2).		OECD recommend that food rations
			increase from 30 mg/larva/day at test
			initiation (Study Day 0-4) to 80
			mg/larva/day in the last week of the test
			(Study Day 15-21).

Page 11 of 61

DER Template Version: 22 September 2011

# Data Evaluation Record on the Toxicity of 2,4-D to Amphibians, Metamorphosis Assay

## EPA MRID Number 48317002

#### C. Exposure System

Table 4: Summary of Information on the Exposure System and Test Vessel Characteristics.

Parameter	Value(s)	Details or Remarks	Guideline Recommendations
Type of exposure:	flow-through		EPA recommends the use of a flow-through system.
Type of flow-through dilution system (if applicable):	continuous-flow diluter system		Intermittent flow proportional diluters or continuous flow serial diluters are recommended.²
Flow-through rate (if applicable):	32 mL/min	± 3 (complete volume turnover approximately every 2.7 hours for each replicate test vessel)	Recommended flow-through rate is 25 mL/min (complete volume replacement ca. every 2.7 hrs).
Details on toxicant mixing for flow-through systems (if applicable):	Stainless steel mixing chamber was used.	The dilution water flow and resulting test solution volume in the mixing chamber was maintained using a stainless steel float valve. The test solution was gravity fed	Recommended toxicant mixing for flow-through systems: 1) Mixing chamber is recommended but not required; 2) Aeration is not recommended for mixing;

Additional guidance for aquatic test design is located in OCSPP Guideline 850.1000, Special Considerations for Conducting Aquatic Laboratory Studies.

Page 12 of 61

DER Template Version: 22 September 2011

Data Evaluation Record on the Toxicity of 2,4-D to Amphibians, Metamorphosis Assay

Parameter	Value(s)	Details or Remarks	Guideline Recommendations
		from the mixing chamber and split equally to replicate test vessels (contained in a temperature controlled water trough) via a Teflon manifold with Teflon delivery tubing. The diluter system was calibrated prior to test initiation, and diluter operation was monitored at least twice daily throughout the test. If there was an indication that the diluter calibration had changed (e.g., diluter malfunction, change in water quality or measured concentrations), calibration of necessary diluter components was checked and adjusted as appropriate.	3) A demonstration that the test solution is completely mixed before introduced into the test system is recommended; 4) The recommended flow splitting accuracy is within 10%.
Renewal period for static renewal (if applicable):	NA		If static renewal is used, EPA recommends 24-hr renewal; renewal period is recommended not to exceed 72 hours.

Page 13 of 61

DER Template Version: 22 September 2011

Data Evaluation Record on the Toxicity of 2,4-D to Amphibians, Metamorphosis Assay

Parameter	Value(s)	Details or Remarks	Guideline Recommendations
Aeration?	No		EPA recommends maintaining dissolved
			oxygen concentrations >40% air
			saturation (>3.5 mg/L). Aeration may
			be maintained through bubblers. It is
			recommended to set bubblers at levels
			that do not cause stress on the tadpoles.
Source of dilution water:	natural water	Water was obtained from the upper	EPA recommends natural or reconstituted
		Saginaw Bay of Lake Huron off Whitestone	water; it is recommended that natural
		Point by the City of Midland Water	water be sterilized with UV and tested for
		Treatment Plant and supplied to The Dow	pesticides, heavy metals, and other
		Chemical Company prior to municipal	possible contaminants, including known
		treatment for human consumption. Water	substrates of the iodine transporter of the
		was limed and flocculated with ferric	thyroid gland (e.g., fluoride, chlorate,
		chloride, sand-filtered, pH-adjusted with	perchlorate). OECD accepts any water in
		gaseous CO2, carbon-filtered, and UV-	which the test species show control
		irradiated.	survival at least as good as indicated in
			the test guideline.

Page 14 of 61

DER Template Version: 22 September 2011

Data Evaluation Record on the Toxicity of 2,4-D to Amphibians, Metamorphosis Assay

Parameter	Value(s)	Details or Remarks	Guideline Recommendations
Was dilution water analyzed for pesticides, heavy metals, and other contaminants?	Yes	Laboratory dilution water is routinely analyzed for pesticides, organics, metals, and other inorganics twice per year. In particular, chlorate, perchlorate, and iodide levels are also determined.	
lodide supplementation in water?	OZ		If reconstituted water is used or if background levels of iodide in natural water are less than 0.5 µg/L, iodide supplementation is recommended. This supplementation is in addition to the recommended dietary source of iodide (e.g. in Sera Micron).
Test vessel type / materials:	Glass sealed together with clear silicone adhesive		EPA and OECD recommend that water-contact portions of the system not compromise the study (e.g., all glass vessels or glass vessels with stainless steel frames are acceptable examples).

Page 15 of 61

DER Template Version: 22 September 2011

Data Evaluation Record on the Toxicity of 2,4-D to Amphibians, Metamorphosis Assay

Parameter	Value(s)	Details or Remarks	Guideline Recommendations
Test vessel size:	30 x 14.5 x 20 cm deep		
Fill volume:	5.2 L		
Additional details on exposure			

Table 5: Summary of Water Quality Characteristics in the Test System.

Parameter	Minimum	Maximum	Mean	Measurement Interval	Guideline Recommendations
Hardness (mg/L as CaCO <sub>3</sub> )	64	92	29	7 days	EPA recommends hardness 40 to 48 mg/L as $CaCO_3$ .
Hd	7.0	7.8		7 days	EPA recommends pH 7.5 $\pm$ 1, interreplicate and inter-treatment differentials should not exceed 0.5.

Page 16 of 61

DER Template Version: 22 September 2011

Data Evaluation Record on the Toxicity of 2,4-D to Amphibians, Metamorphosis Assay

Parameter	Minimum	Maximum	Mean	Measurement	Guideline Recommendations
Dissolved oxygen (mg/L)	5.3	7.9		7 days	EPA recommends dissolved oxygen (DO) >3.5 mg/L (>40% air saturation). OECD recommends DO concentration >3.5 mg/L (>40% air saturation).
Temperature	21.9	22.5		7 days	EPA recommends temperature $22\pm 1^{\circ}C$ ; inter-replicate and inter-treatment differentials should not exceed $0.5^{\circ}C$ .
lodide	CL0Q	<07>			EPA recommends aquatic iodide range 0.5 - 10 µg/L (supplemental iodide should not exceed 2 µg/L).
Ammonia	<0.1	0.24		4 months	General recommendations for frequency
Fluoride	0.1	0.1			of measurements: EPA recommends that
Perchlorate	<0.2	<0.2			water quality parameters be measured in a control and at one test item
Chlorate	<10	<10			concentration at least weekly. In static renewal systems, water quality

Page 17 of 61

DER Template Version: 22 September 2011

Data Evaluation Record on the Toxicity of 2,4-D to Amphibians, Metamorphosis Assay

Parameter	Minimum	Maximum	Mean	Measurement Interval	Guideline Recommendations
					parameters, including ammonia, should
					be measured just prior to renewal. In
					addition, EPA recommends that DO be
					measured at each concentration at least
					weekly and that temperature be
					measured continuously. OECD
					recommends that DO and temperature
					be measured at least weekly and that pH
					and hardness be measured at least at
					the beginning and end of the test.

Page 18 of 61

DER Template Version: 22 September 2011

D. Study Design and Additional Experimental Conditions

Table 6: Range-Finding Study Conditions (if Applicable).

Parameter	Value(s)	Details or Remarks	Guideline Recommendations
Was a range-finder conducted?	ON		
If yes, what was the method for	AN		EPA recommends that the highest test
determining the highest test			concentration is either the solubility limit of
concentration in the range-			the test compound, 100 mg/L, or
finder?			demonstrates adequate evidence of toxicity
			(e.g., <10% mortality), whichever
			concentration is lowest.
Species:	NA		
Life stage:	NA		
Test duration:	NA		
Additional details:	NA		

Page 19 of 61

DER Template Version: 22 September 2011

# Data Evaluation Record on the Toxicity of 2,4-D to Amphibians, Metamorphosis Assay

Table 7: Definitive Study Conditions.

Parameter	Value(s)	Details or Remarks	Guideline Recommendations
Test duration:	21 days		EPA recommends that the duration of the definitive test be 21 days.
Method for selecting the highest test concentration in the definitive test:	The solubility limit and other reference studies were used to select the highest test concentration.		EPA recommends that the highest test concentration is either the solubility limit of the test compound, 100 mg/L, or demonstrates adequate evidence of toxicity (e.g., <10% mortality), whichever concentration is lowest.
Reference study citation (if applicable):	Morgan et al., 1996; Palmer and Krueger, 1997a; Palmer and Kreuger, 1997b; Alexander et al., 1985		
Separation of test concentrations:	0.1-0.4		EPA recommends that the maximum concentration separation be 0.1 and the minimum be 0.33.

Page 20 of 61

DER Template Version: 22 September 2011

Data Evaluation Record on the Toxicity of 2,4-D to Amphibians, Metamorphosis Assay

Parameter	Value(s)	Details or Remarks	Guideline Recommendations
Number of test concentrations:	4		EPA recommends a minimum of 3 concentrations and a control, plus solvent control if appropriate.
Are nominal concentrations adjusted for purity?	O <sub>Z</sub>	"Standards were not corrected for purity." (Appendix B)	
Indicate the type of values presented for measured concentrations:	Average measured		
Limit of quantification (LOQ):	0.120 mg a.i./L		EPA recommends that for chemical test concentrations below the LOQ, analyses be conducted on the stock solutions.
Level of detection (LOD):	Not reported		
Frequency of measurement:	7 days		It is recommended that test item concentration be measured in one tank at each treatment level at test initiation and every week thereafter.

Page 21 of 61

DER Template Version: 22 September 2011

Data Evaluation Record on the Toxicity of 2,4-D to Amphibians, Metamorphosis Assay

Parameter Value(s)	Value(s)	Details or Remarks	Guideline Recommendations
Number of replicates in control:	4		EPA recommends 4 replicates.
Number of replicates in solvent			EPA and OECD recommend the use of a
control (if applicable):			concurrent solvent control when a
			solubilizing agent is used. EPA
			recommends 4 replicates.
Number of replicates per test	4		EPA recommends 4 replicates.
item treatment level:			
Number of larvae per treatment	80		
at test initiation:			
Was a solvent used?	o Z		
Solvent type (if applicable):			
Maximum solvent concentration	Ϋ́		EPA recommends that the solvent not
(if applicable):			exceed 0.02 ml/L³. OECD recommends

Hutchinson TH, Shillabeer N, Winter MJ, Pickford DB (2006). Acute and chronic effects of carrier solvents in aquatic organisms: A critical review. Review. Aquatic Toxicology, 76, pp.69-92.

Page 22 of 61

DER Template Version: 22 September 2011

Data Evaluation Record on the Toxicity of 2,4-D to Amphibians, Metamorphosis Assay

Parameter	Value(s)	Details or Remarks	Guideline Recommendations
			that solvent have no effect on survival nor produce any other adverse effects and that
Was a positive control used?	ON		concentration not be greater than U.1 mi/L.
Positive control (if applicable):	AN		
Positive control concentration(s) (if applicable):	NA		
Photoperiod:	12 hrs light : 12 hrs dark		EPA recommends photoperiod 12:12 (light:dark).
Light intensity at water's surface:	0.618-0.882 Klux		EPA recommends light intensity 0.6 - 2 Klux (at water's surface).

OECD (2000). Guidance Document on Aquatic Toxicity Testing of Difficult Substances and Mixtures. Environmental Health and Safety Publications. Series on Testing and Assessment. No. 23. Paris, France.

Page 23 of 61

DER Template Version: 22 September 2011

Data Evaluation Record on the Toxicity of 2,4-D to Amphibians, Metamorphosis Assay

Parameter	Value(s)	Details or Remarks	Guideline Recommendations
Additional details:		Information on test solution appearance	
		did not appear to be reported in the	
		report.	

Table 8: Summary of Treatment Concentrations in the Amphibian Metamorphosis Assay with 2,4-D.

Guideline Recommendations	EPA and OECD recommend that test item	concentrations be maintained at a coefficient	or variation (CV) SCO%.		
Details or Remarks					
Mean CV (%)		55.67	21.60	16.29	16.90
Measured Concentration (mg a.i./L)	<000	0.273	3.24	38.0	113
Nominal Concentration (mg a.i./L)	0	0.4	4	40	100
Treatment ID	Negative Control	Treatment 1	Treatment 2	Treatment 3	Treatment 4

Abbreviations: CV Coefficient of variation.

LOQ = 0.120 mg a.i./L

Page 24 of 61

DER Template Version: 22 September 2011

# Data Evaluation Record on the Toxicity of 2,4-D to Amphibians, Metamorphosis Assay

#### EPA MRID Number 48317002

#### E. Observations

Mortality (throughout study); NF stage and asynchronous development; wet weight, snout-vent length, hind limb length (Day Biological Endpoints:

7 and 21); thyroid histopathology (Day 21)

Were raw (individual) data provided? Yes

EPA recommends that observations of mortality and clinical signs occur daily, at a minimum; other observations are recommended as follows: NF developmental stage (Days 7 and 21); any asynchronous development, indicated by tadpoles that cannot be assigned an NF stage (Days 7 and 21); hind limb length (Days 7 and 21); snout-vent length (Days 7 and 21); body weight (test initiation, for optional size-based larval selection); and thyroid gland gross pathology and histopathology (Day 21). Note the histopathology section of the test guideline also includes thyroid gross pathology observations.

Page 25 of 61

DER Template Version: 22 September 2011

#### II. RESULTS AND DISCUSSION

#### A. Results

Only one larvae in the 3.24 mg a.i./L treatment level died during the 21-day study.

Table 9: Larval Mortality in South African Clawed Frog.

			Larval N	/lortality	,	
Treatment (mg a.i./L) [mean-measured]		Day 7 <sup>1</sup>			Day 21	
[moun moudulou]	n	Mortality #	Mortality %	n	Mortality #	Mortality %
Negative Control	80	0	0	60	0	0
0.273	80	0	0	60	0	0
3.24	80	0	0	60	1	1.7
38.0	80	0	0	60	0	0
113	80	0	0	60	0	0

Sample size and cumulative mortality values at Day 7 prior to interim sacrifice.

Day 7 median NF stage was 54 for the negative control and all treatment levels. Day 21 median NF stage was 58 for the highest treatment level and 59 for the negative control and all other treatment levels. There were no asynchronous tadpoles.

Table 10: Larval Development in South African Clawed Frog- Developmental Stage and Asynchronous Development.

			Developme	ental St	age	
Treatment (mg a.i./L)		Da	ny 7		Da	y 21
[mean-measured]	n	Median Stage <sup>1</sup>	# Asynchronous	n	Median Stage <sup>1</sup>	# Asynchronous
Negative Control	4	54	0	4	59	0
0.273	4	54	0	4	59	0
3.24	4	54	0	4	59	0
38.0	4	54	0	4	59	0
113	4	54	0	4	58	0

Page 27 of 61

DER Template Version: 22 September 2011

Day 7 HLL ranged from 2.2 to 2.4 mm across the control and all treatment levels. Day 21 HLL ranged from 16.2 mm in the high treatment group to 18.0 mm in the low treatment group.

Table 11: Larval Development in South African Clawed Frog - Hind Limb Length.

				Hind Limb	Lengt	h (HLL)		
Treatment			Day 7				Day 21	
(mg a.i./L) [mean-measured]	n	Mean (mm)	±SD	HLL: SVL <sup>1</sup>	n	Mean (mm)	±SD	HLL: SVL <sup>1</sup>
Negative Control	4	2.4	0.2	0.1	4	17.2	0.8	0.7
0.273	4	2.2	0.1	0.1	4	18.0	1.0	0.7
3.24	4	2.3	0.3	0.1	4	17.1	0.9	0.7
38.0	4	2.2	0.1	0.1	4	17.5	1.0	0.7
113	4	2.2	0.1	0.1	4	16.2	1.0	0.6

Abbreviations: SD Standard deviation.

In this table, "n" represents the number of independent replicates per treatment level.

Summary results for snout-vent length (SVL) are presented in the next table (Table 12).

Day 7 SVL ranged from 18.00 to 18.5 mm and Day 21 SVL ranged from 26.9 to 28.1 mmm across the negative control and all treatment levels. Day 7 body weight ranged from 0.48 g to 0.53 g and Day 21 body weight ranged from 1.63 g to 1.91 g across the negative control and all treatment levels.

Table 12: Larval Growth in South African Clawed Frog.

Tuestusent		Sno	ut-Vent	Length	ı (SVL)				Body V	Veigl	nt¹	
Treatment (mg a.i./L)		Day 7			Day 21			Day 7	7		Day 2	11
[measured]	n	Mean (mm)	±SD	n	Mean (mm)	±SD	n	Mean (g)	±SD	n	Mean (g)	±SD
Negative Control	4	18.5	0.5	4	26.8	1.4	4	0.53	0.04	4	1.63	0.17
0.273	4	18.0	0.5	4	27.4	0.6	4	0.48	0.04	4	1.73	0.12
3.24	4	18.3	0.6	4	26.9	1.6	4	0.49	0.05	4	1.63	0.24
38.0	4	18.1	0.2	4	27.7	0.31	4	0.49	0.02	4	1.75	0.06
113	4	18.4	0.3	4	28.1	0.40	4	0.50	0.04	4	1.91	0.22

Abbreviations: SD Standard deviation.

In this table, "n" represents the number of independent replicates per treatment level.

There were no treatment-related histopathologic changes in the thyroid gland in any of the treatment groups. There was no evidence of glandular atrophy or hypertrophy or follicular cell hyperplasia in any of the thyroid glands examined across all treatment groups. The incidence of tall columnar cells lining the follicles (follicular cell hypertrophy) did not show any treatment-related differences and was interpreted to be within normal limits at all concentrations of 2,4-D.

Also referred to as "wet weight" in the test guideline.

Table 13: Gross Pathology and Histopathology of the Thyroid Gland in South African Clawed Frog.

Treatment				Diagr	nostic Obser	vations	<b>S</b> <sup>1</sup>		
(mg a.i./L) [mean-	Severity	-	oid Gland pertrophy		roid Gland Atrophy		icular Cell pertrophy		cular Cell erplasia
measured]		n	Incidence	n	Incidence	n	Incidence	n	Incidence
Negative	0	20	20	20	20	20	5	20	20
Control	1	20	0	20	0	20	15	20	0
	2	20	0	20	0	20	0	20	0
	3	20	0	20	0	20	0	20	0
0.273	0	20	20	20	20	20	3	20	20
	1	20	0	20	0	20	17	20	0
	2	20	0	20	0	20	0	20	0
	3	20	0	20	0	20	0	20	0
3.24	0	20	20	20	20	20	6	20	20
	1	20	0	20	0	20	14	20	0
	2	20	0	20	0	20	0	20	0
	3	20	0	20	0	20	0	20	0
38.0	0	20	20	20	20	20	5	20	20
	1	20	0	20	0	20	15	20	0
	2	20	0	20	0	20	0	20	0
	3	20	0	20	0	20	0	20	0
113	0	20	20	20	20	20	4	20	20
	1	20	0	20	0	20	16	20	0
	2	20	0	20	0	20	0	20	0
	3	20	0	20	0	20	0	20	0

Thyroid gland gross pathology and histopathology are graded 0 - 3 based on severity: 0=Not remarkable, 1=Mild, 2=Moderate, 3=Severe. See OECD No. 82 for reference.

Page 30 of 61

DER Template Version: 22 September 2011

Data Evaluation Record on the Toxicity of 2,4-D to Amphibians, Metamorphosis Assay

Table 14: Additional Thyroid Gland Histopathology Observations in South African Clawed Frog.

Treatment					Additional Qu	ıalitative	Additional Qualitative Observations				
(mg a.i./L)		Follic	Follicular Lumen	Folli	Follicular Lumen	Foll	Follicular Cell	Follicul	Follicular Cell Height	Fol	Follicular Cell
[time-weighted,	Severity	Area	Area (Increase)	Area	Area (Decrease)	Heigł	Height (Increase)	Q)	(Decrease)		Shape
mean-measured]		u	Incidence	п	Incidence	п	Incidence	u	Incidence	u	Incidence
Negative Control	0	ΑN	AN	Ą	AN	AN	ΑN	NA	ΑN	Ϋ́Z	NA
	1	ΝA	ΥN	NA	VΝ	NA	NA	NA	AN	NA	NA
	2	ΝA	ΥN	NA	VΝ	NA	NA	NA	AN	NA	NA
	3	NA	NA	NA	NA	NA	NA	NA	NA	NA	NA
0.273	0	NA	NA	NA	NA	NA	NA	NA	NA	NA	NA
	1	ΝA	ΝA	NA	VΝ	NA	NA	NA	NA	NA	NA
	2	NA	NA	NA	NA	NA	NA	NA	NA	NA	NA
	3	NA	NA	NA	NA	NA	NA	NA	NA	NA	NA
3.24	0	NA	NA	NA	NA	NA	NA	NA	NA	NA	NA
	1	NA	NA	NA	NA	NA	NA	NA	NA	NA	NA
	2	NA	NA	NA	NA	NA	NA	NA	NA	NA	NA
	3	NA	NA	NA	NA	NA	NA	NA	NA	NA	NA

Page 31 of 61

DER Template Version: 22 September 2011

Data Evaluation Record on the Toxicity of 2,4-D to Amphibians, Metamorphosis Assay

EPA MRID Number 48317002

Treatment					Additional Qu	ıalitative	Additional Qualitative Observations				
(mg a.i./L)		Follic	Follicular Lumen	Folli	Follicular Lumen	Fol	Follicular Cell	Follicul	Follicular Cell Height	Fol	Follicular Cell
[time-weighted,	Severity	Area	Area (Increase)	Area	Area (Decrease)	Heig	Height (Increase)	<u>O</u>	(Decrease)		Shape
mean-measured]		и	Incidence	u	Incidence	u	Incidence	и	Incidence	u	Incidence
38.0	0	NA	NA	NA	NA	NA	NA	NA	NA	NA	NA
	1	NA	NA	NA	NA	NA	NA	NA	NA	NA	NA
	2	NA	NA	AN	NA	NA	NA	NA	NA	NA	NA
	3	NA	NA	NA	NA	NA	NA	NA	NA	NA	NA
113	0	NA	NA	NA	NA	NA	NA	NA	NA	NA	NA
	1	NA	NA	NA	NA	NA	NA	NA	NA	NA	NA
	2	NA	NA	NA	NA	NA	NA	NA	NA	NA	NA
	3	NA	NA	NA	NA	NA	NA	NA	NA	NA	NA

Abbreviations: Not applicable.

Page 32 of 61

DER Template Version: 22 September 2011

Thyroid histopathology is graded 0 - 3 based on severity: 0=Not remarkable, 1=Mild, 2=Moderate, 3=Severe. See OECD No. 82 for reference.

No abnormal behaviors or clinical signs of toxicity were noted among control or 2,4-D exposed tadpoles.

Table 15: Clinical Signs in Xenopus laevis.

Treatment	Clinical Signs <sup>1</sup>		
(mg a.i./L)	Туре	n	Incidence
[TWA- measured]	туре	"	incidence
Negative Control	None	NA	NA
Solvent Control	None	NA	NA
0.00264	None	NA	NA
0.0262	None	NA	NA
0.273	None	NA	NA

Abbreviations: NA Not applicable.

Note that asynchronous development (unable to stage) is reported previously in Table 10 and not here.

B. Study Author's Analysis and Conclusions

There were no signs of overt toxicity among exposed tadpoles in the present study. Throughout the

entire exposure period, there was only one incidence of tadpole mortality. The study author identified

no indications of developmental delay or abnormal behavior and therefore concluded that concentrations

of 2,4-D used in the present study (up to 113 mg/L 2,4-D) were not overtly toxic to developing X.

laevis. There were also no signs of advanced development (as measured by developmental stage and

hind limb length) or asynchronous development among 2,4-D exposed tadpoles relative to control

tadpoles on either Day 7 or Day 21 of exposure. Finally, compared to thyroid glands from controls,

there were no significant histopathological effects observed among thyroid glands from 2,4-D exposed

tadpoles.

C. Reviewer's Analysis and Conclusions

Statistical Methods: Day 21 values for NF developmental stage and normalized HLL were consistent

with a monotonic response and were analyzed using the Jonckheere-Terpstra test. Values for Day

7 SVL and Days 7 and 21 body weight showed no apparent monotonic response and met the

assumptions of normality and homogeneity of variance; therefore, these endpoints were analyzed using

the parametric Dunnett's test. Day 7 NF stage and Day 21 SVL were inconsistent with a monotonic

response and were analyzed with the non-parametric Mann-Whitney test. All analyses were performed

using the FROG program with Statistical Analysis Software (SAS v. 9.3).

Late stage (> NF 60) tadpoles were excluded from analyses of SVL, body weight, and normalized

HLL. Late stage tadpoles were excluded from every treatment level: control-6 tadpoles, 0.273 mg

a.i./L-10 tadpoles, 3.25 mg a.i./L-6 tadpoles, 38.0 mg a.i./L-7 tadpoles, 113 mg a.i./L-6 tadpoles.

Histopathology results were evaluated visually based on severity and incidence data, within the context

of the narrative pathology report.

Page 34 of 61

DER Template Version: 22 September 2011

# Conclusions:

The only statistically-significant effect was for Day 21 HLL which was decreased by 15% at the high treatment level compared to the negative control after excluding "late stage" tadpoles from the analysis.

Page 35 of 61

# Data Evaluation Record on the Toxicity of 2,4-D to Amphibians, Metamorphosis Assay

# EPA MRID Number 48317002

Table 16: Developmental and Thyroid Gross Pathology/Histopathology Endpoints<sup>1,2</sup> in the AMA with 2,4-D.

	Ц	NE Develonmental	ental Stade		I	mi I bui	Hind Limb Length <sup>3</sup>	3		Asynch	Asynchronous		Thyroid Gross and
Treatment										Develo	Development		Histopathology
(mg a.i./L)	Day 7	7	Day 21	21	Day 7	, 7	Day 21	21	Day 7	7 /	Day 21	. 21	Day 21
[mean-measured]	: ( N	1	: ( ) N	1	%	1	%	1	%	1	%	1	Treatment-Related
	Median	р	Median	р	Diff.	р	Diff.	р	Diff.	р	Diff.	ď	Effects? (Yes/No)
Negative Control	54	NA	65	NA	NA	NA	NA	NA	NA	NA	NA	NA	NA
0.273	54	0.48	65	0.75	0	AN	-3.7	0.25	0	ΑN	0	A N	ON.
3.24	54	>0.99	65	0.94	0	NA	-3.7	0.25	0	NA	0	NA	No
38.0	54	>0.99	65	0.59	0	NA	-3.7	0.26	0	NA	0	NA	No
113	54	>0.99	28	0.17	0	NA	-14.8	0.02	0	NA	0	NA	No
Statiotical Toot	Month distribution	, o atial	Jonckheere-	ere-	Ž		Jonckheere-	leere-	VIZ	<	2	<u> </u>	Ž
Statistical Test	Malli	niii iey	Terpstra	tra	2	ſ	Terpstra	stra	2	(	2	ζ	<u> </u>

Abbreviations: Difference. NA Not applicable.

Unless otherwise indicated, effects are reported based on comparison to the clean water control. Negative values indicate reductions or decreases relative to the negative control. Conclusions regarding histopathology may be heavily weighted by the expert opinion of a board-certified pathologist.

<sup>2</sup> Unless otherwise specified, effects are considered statistically significant at p<0.05.</p>

Hind-limb length is normalized to snout-vent length (SVL).

Page 36 of 61

Table 17: Growth Endpoints<sup>1,2</sup> in the AMA with 2,4-D.

Tuestment		Snout-Ve	nt Length			Body V	Veight	
Treatment	Day	7	Day	21	Day	7	Day	21
(mg a.i./L) [mean-measured]	% Diff.	р	% Diff.	р	% Diff.	р	% Diff.	р
Negative Control	NA	NA	NA	NA	NA	NA	NA	NA
0.273	-2.8	0.31	2.3	0.78	-9.0	0.32	6.3	0.82
3.24	-1.1	0.92	0.3	>0.99	-6.3	0.61	-0.8	>0.99
38.0	-2.2	0.53	3.4	0.78	-7.2	0.50	7.9	0.68
113	-0.7	0.98	4.9	0.24	-5.8	0.68	17.5	0.10
Statistical Test	Dunn	ett's	Mann-V	Vhitney	Dunne	ett's	Dunr	nett's

Abbreviations: Difference. NA Not applicable.

# E. Study Deficiencies

The coefficients of variation for the measured concentrations of the 0.273 and 3.24 mg a.i./L treatment groups exceeded the guideline performance criterion of 20%. This was likely due to biodegradation of the test material in the test vessels. All other validity and performance criteria were met.

# F. Reviewer's Comments

Unlike the study author, the reviewer detected a statistically significant decrease in Day 21 normalized HLL (p<0.05) at the highest treatment level after excluding "late stage" tadpoles from the analysis. No other signs of toxicity were observed at this or any other treatment level. Late stage tadpoles were excluded from every treatment level: control-6 tadpoles, 0.273 mg a.i./L-10 tadpoles, 3.25 mg a.i./L-6 tadpoles, 38.0 mg a.i./L-7 tadpoles, 113 mg a.i./L-6 tadpoles. Late stage tadpoles were included in the analysis of potential effects on developmental stage.

Page 37 of 61

Unless otherwise indicated, effects are reported based on comparison to the negative (clean water) control.

<sup>&</sup>lt;sup>2</sup> Unless otherwise specified, effects are considered statistically significant at p<0.05.

The individuals selected for histopathology were not stage-matched to the median stage in the controls. However, there were no noteworthy differences between histopathology observations in control and treatment specimens; therefore, this guideline deviation does not appear to have substantively affected the interpretation of results.

### III. REFERENCES

Alexander, H. C., Gersich, F. M., and Mayes, M. A. (1985). Acute Toxicity of Four Phenoxy Herbicides to Aquatic Organisms. Bulletin of Environmental Contamination and Toxicology 35, 314-321.

ASTM. (2004). Standard Guide for Conducting the Frog Embryo Teratogenesis Assay - Xenopus (FETAX). E 1439-98

ASTM. (2007). Standard Guide for Conducting Acute Toxicity Tests on Test Materials with Fishes, Macroinvertebrates, and Amphibians. American Society for Testing and Materials, ASTM £729-96(2002), Philadelphia, Pennsylvania.

Cooke, A. S. (1972). The effects of DDT, Dieldrin, and 2,4-D on amphibian spawn and tadpoles. Environ. Pollut. 3, 51-68.

ECOTOX database. 2009. Search of the United States Environmental Protection Agency's ECOTOX database on CAS# 94757 and Pimephales promelas conducted on 12/17/09; http://cfpub.epa.gov/ecotox/.

Hansch C, Leo A, and Hoekman D. (1995). Exploring QSAR Hydrophobic, Electronic and Steric Constants, ACS, Washington, D.C.

Hollander, M. and Wolfe, D. A. (1973). Nonparametric Statistical Methods. John Wiley & Sons, New York.

Page 38 of 61

- Hopkins, D. L. (1987). 2,4-Dichlorophenoxyacetic Acid: Determination of the water solubility. Unpublished report of The Dow Chemical Company. Laboratory Report No: ES-994. November 25, 1987. MRID No. 41332002.
- Lezotte F., Van Hoven R., and Nixon, W. (2002). Determination of Water Solubility of 2,4- D Acid by the Shake Flask Method: Lab Project Number: 467C-103: 467/121801/7840/SUB467. Unpublished study prepared by Wildlife International, Ltd. MRID No. 45692501.
- Morgan, M. K., Scheuerman, P. R., Bishop, C. R., and Pyles, R. A. (1996). Teratogenic potential of atrazine and 2,4-D using FETAX. Journal of Toxicology and Environmental Health 48, 151-168.
- Nieuwkoop, P. D. and Faber, J. (1994). Normal Table of Xenopus laevis (Daudin). Garland Publishing, New York
- Organization for Economic Cooperation and Development (OECD). (2000). Guidance Document on Aquatic Toxicity Testing of Difficult Substances and Mixtures. Environmental Health and Safety Publications. Series on Testing and Assessment. No. 23. Paris.
- Organization for Economic Cooperation and Development OECD. (2007). Guidance Document on Amphibian Thyroid Histology. Environmental Health and Safety Publications. Series on Testing and Assessment. No. 82. Paris
- Organization for Economic Cooperation and Development (OECD). (2009). OECD Guideline for the Testing of Chemicals: The Amphibian Metamorphosis Assay. OECD 231 Adopted 7 September 2009.
- Palmer, S. J. and Krueger, H. O. (1997a). 2,4-D Dimethylamine salt: A 96-hour static acute toxicity test with the leopard frog tadpoles (Ranapipiens). Study ID 467-103. Wildlife International Ltd., Easton, Maryland. MRID No. 44517306.

Page 39 of 61

- Palmer, S. J. and Krueger, H. O. (1997b). 2,4-D (2,4-Dichlorophenoxyacetic acid): A 96-hour static acute toxicity test with the leopard frog tadpoles (Ranapipiens). Study ID 467-102. Wildlife International Ltd., Easton, Maryland. MRID No. 44517307.
- Sinton, G. L., Fan, L. T., Erickson, L. E., and Lee, S. M. (1986). Biodegradation of 2,4-D and related xenobiotic compounds. Enzyme Microb. Technol. 8, 395-403.
- United States Environmental Protection Agency (USEPA) (2005). Reregistration Eligibility Decision for 2,4-D. EPA 738-R-05-002, June 2005.
- United States Environmental Protection Agency (USEPA) (2009). Endocrine Disrupter Screening Program Test Guidelines OPPTS 890.1100: Amphibian Metamorphosis (Frog). EPA 740-C-09-002, October 2009.
- Wo C. (2008). Certificate of Analysis for 2,4-D Acid, Lot# 2006 2433 8006-USA. PSL Reference No. 081210-2H. Eurofins [Product Safety Laboratories, Dayton, New Jersey.

# Appendix I: Output of Reviewer's Statistical Analysis

test for amphib metamorph screen study - 2,4-d ANALYSIS RESULTS FOR VARIABLE VAR01 ( 7-d wet weight (g) )

TESTS OF ASSUMPTIONS FOR PARAMETRIC ANALYSIS
Shapiro-Wilks test for Normality of Residuals -- alpha-level=0.01
Levenes test for homogeneity of variance(absolute residuals) -- alpha-level=0.05

Levenes test for homogeneity of variance(absolute residuals) -- alpha-level=0.0 Use parametric analyses if neither test rejected, otherwise non-parametric analyses.

Shapiro-Wilks Shapiro-Wilks Levenes Levenes Conclusion

Snapr.	T.O-MT	irks Snat	DITO-WIIKS	Levenes	Levenes	Concrus.	LOH	
Tes	t Sta	at P-	-value	Test Sta	t P-value			
0	.909	(	0.061	0.655	0.632	USE PARA	AMETRIC	TESTS
*****	****	*******	*****	*****	*****	*****	*****	*****
BASIC S	UMMAR	RY STATIST	TICS					
Level	N	Mean	StdDev	StdErr	Coef of Va	r 95%	Conf.In	nterval
Ctrl	4	0.53	0.04	0.02	8.03	(	0.46,	0.59
Dose1	4	0.48	0.04	0.02	8.86	(	0.41,	0.55
Dose2	4	0.49	0.05	0.02	9.97	(	0.41,	0.57
Dose3	4	0.49	0.02	0.01	4.57	(	0.45,	0.52
Dose4	4	0.50	0.04	0.02	8.41	(	0.43,	0.56
Level		Median	Min	Max	%of Control(m	eans) <sup>9</sup>	Reduct:	ion(means)
Ctrl		0.52	0.48	0.58				
Dose1		0.47	0.44	0.54	90.97		9.03	3
Dose2		0.47	0.46	0.57	93.68		6.32	2
Dose3		0.49	0.46	0.51	92.78		7.22	2
Dose4		0.49	0.46	0.54	94.25		5.7	5

\*

PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests Analysis of Variance (ANOVA) - overall F-test

Numerator df Denominator df F-stat P-value 4 15 0.78 0.556

Dunnett - testing each trt mean signif. different than control Williams - test assumes dose-response relationship, testing negative trend Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC

Tues	
Dose4 Do	se5
	•

\*

NON-PARAMETRIC ANALYSES – use alpha-level=0.05 for all tests Kruskal-Wallis test – equality among treatment groups

Page 41 of 61

Degrees of Freedom TestStat P-value 3.62 0.460 MannWhit - testing each trt median signif. different from control Jonckheere - test assumes dose-response relationship, testing negative trend Level Median MannWhit p-value Jonckheere p-value 0.52 Ctrl Dose1 0.47 0.235 0.074 0.47 Dose 2 0.233 0.120 Dose3 0.49 0.230 0.200 0.49 0.494 0.275 Dose4 LOWEST CONCENTRATION SIGNIF. LESS THAN CONTROL DECREASING TREND TEST SUMMARY >highest dose (no sign. differences) Williams Jonckheere >highest dose (no sign. differences) \* PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests Analysis of Variance (ANOVA) - overall F-test Numerator df Denominator df F-stat 15 0.78 0.556 4 Dunnett - testing each trt mean signif. different than control Williams - test assumes dose-response relationship, testing INCREASING trend Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC Dunnett Isotonic Williams Level Tukey p-values mean p-value Dose1 Dose2 Dose3 Dose4 p-value Dose5 -0.53 -0.50 Ctrl -0.50 0.906 -0.50 0.927 -0.48 0.320 Dose1 0.986 Dose2 -0.49 0.611 -0.49 0.504 -0.50 0.937 0.997 1.000 Dose3 Dose4 -0.50 0.680 -0.50 0.943 0.973 1.000 \* NON-PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests Kruskal-Wallis test - equality among treatment groups Degrees of Freedom TestStat P-value 3.62 0.460 MannWhit - testing each trt median signif. different from control Jonckheere - test assumes dose-response relationship, testing INCREASING trend MannWhit p-value Level Median Jonckheere p-value Ctrl -0.52 0.235 0.926 Dose1 -0.47 Dose2 -0.47 0.233 0.880 0.800 Dose3 -0.49 0.230 Dose4 -0.49 0.494 0.725

Page 42 of 61

```
INCREASING TREND TEST SUMMARY LOWEST CONCENTRATION SIGNIF. GREATER THAN
   Williams
                                                          >highest dose (no sign. differences)
   Jonckheere
                                                          >highest dose (no sign. differences)
test for amphib metamorph screen study - 2,4-d
ANALYSIS RESULTS FOR VARIABLE VAR02 ( 7-d stage (median) )
TESTS OF ASSUMPTIONS FOR PARAMETRIC ANALYSIS
Shapiro-Wilks test for Normality of Residuals -- alpha-level=0.01
Levenes test for homogeneity of variance(absolute residuals) -- alpha-level=0.05
Use parametric analyses if neither test rejected, otherwise non-parametric
analyses.
  Shapiro-Wilks Shapiro-Wilks
                                         Levenes
                                                       Levenes Conclusion
                                       Test Stat P-value
     Test Stat P-value
                       < .001
                                                        <.001
       0.509
                                         9.000
                                                                    USE NON-PARAMETRIC TESTS
*************************
BASIC SUMMARY STATISTICS
 Level N Mean StdDev StdErr Coef of Var 95% Conf.Interval
  Ctrl 4 54.00 0.00 0.00 0.00

Dosel 4 53.75 0.50 0.25 0.93

Dosel 4 54.00 0.00 0.00 0.00

Dosel 4 54.00 0.00 0.00 0.00

Dosel 4 54.00 0.00 0.00 0.00

Dosel 4 54.00 0.00 0.00 0.00
                                                                          52.95,
                                                                                        54.55

        Median
        Min
        Max
        % of Control(means)
        % Reduction(means)

        54.00
        54.00
        .
        .

        54.00
        53.00
        54.00
        99.54
        0.46

        54.00
        54.00
        100.00
        0.00

        54.00
        54.00
        100.00
        0.00

        54.00
        54.00
        100.00
        0.00

        54.00
        54.00
        100.00
        0.00

           Median
 Level
  Ctrl
  Dose1
  Dose2
  Dose3
  Dose4
******************
PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
     Analysis of Variance (ANOVA) - overall F-test
      Numerator df Denominator df F-stat
                                                                 P-value
                            15
                                                1.00
                                                                0.438
             4
Dunnett - testing each trt mean signif. different than control
Williams - test assumes dose-response relationship, testing negative trend
Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC
Level Mean
                   Dunnett Isotonic Williams
                                                                           Tukey p-values
                                          p-value Dose1 Dose2 Dose3 Dose4
                   p-value
                               mean
                                                                                                 Dose5
          54.00
                                54.00
 Ctrl
 Dosel 53.75 0.357 53.94 0.415
                                                                  •

      Dose2
      54.00
      1.000
      53.94
      0.443
      0.530
      .

      Dose3
      54.00
      1.000
      53.94
      0.458
      0.530
      1.000

      Dose4
      54.00
      1.000
      53.94
      0.468
      0.530
      1.000

                                            0.468 0.530 1.000 1.000
*************************
```

Page 43 of 61

```
NON-PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
    Kruskal-Wallis test - equality among treatment groups
     Degrees of Freedom TestStat P-value
                             4.00
MannWhit - testing each trt median signif. different from control
Jonckheere - test assumes dose-response relationship, testing negative trend
         Median
                          MannWhit p-value
                                                       Jonckheere p-value
 Level
          54.00
  Ctrl
           54.00
                                     0.478
                                                              0.159
  Dose1
           54.00
54.00
                                     1.000
                                                              0.500
  Dose2
  Dose3
                                     1.000
                                                              0.673
           54.00
  Dose4
                                     1.000
                                                              0.760
 DECREASING TREND TEST SUMMARY
                                   LOWEST CONCENTRATION SIGNIF. LESS THAN CONTROL
   Williams
                                               >highest dose (no sign. differences)
   Jonckheere
                                                >highest dose (no sign. differences)
******************
PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
    Analysis of Variance (ANOVA) - overall F-test
     Numerator df Denominator df F-stat
                                                      P-value
                                        1.00
          4
                      15
                                                     0.438
Dunnett - testing each trt mean signif. different than control
Williams - test assumes dose-response relationship, testing INCREASING trend
Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC
       Mean Dunnett Isotonic Williams
Level
                                                              Tukey p-values
                         mean p-value Dose1 Dose2 Dose3 Dose4
                                                                               Dose5
               p-value
 Ctrl
       -54.00
                          -53.88

      Ctrl
      -54.00
      .
      -53.88
      .
      .
      .
      .

      Dosel
      -53.75
      0.357
      -53.88
      0.855
      .
      .
      .

      Dose2
      -54.00
      1.000
      -54.00
      0.618
      0.530
      .

      Dose3
      -54.00
      1.000
      -54.00
      0.636
      0.530
      1.000

 Dose4 -54.00 1.000
                          -54.00
                                    0.648 0.530 1.000 1.000
************************
NON-PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
    Kruskal-Wallis test - equality among treatment groups
     Degrees of Freedom TestStat P-value
                             4.00
                                        0.406
MannWhit - testing each trt median signif. different from control
Jonckheere - test assumes dose-response relationship, testing INCREASING trend
                          MannWhit p-value
 Level
          Median
                                                      Jonckheere p-value
  Ctrl
          -54.00
          -54.00
                                    0.478
                                                             0.841
  Dose1
                                    1.000
                                                             0.500
  Dose2
          -54.00
  Dose3
          -54.00
                                    1.000
                                                             0.327
          -54.00
                                     1.000
                                                              0.240
  Dose4
```

Page 44 of 61

```
INCREASING TREND TEST SUMMARY LOWEST CONCENTRATION SIGNIF. GREATER THAN
CONTROL.
  Williams
                                           >highest dose (no sign. differences)
  Jonckheere
                                          >highest dose (no sign. differences)
test for amphib metamorph screen study - 2,4-d
ANALYSIS RESULTS FOR VARIABLE VAR03 ( 7-d sn-vent length (mm) )
TESTS OF ASSUMPTIONS FOR PARAMETRIC ANALYSIS
Shapiro-Wilks test for Normality of Residuals -- alpha-level=0.01
Levenes test for homogeneity of variance(absolute residuals) -- alpha-level=0.05
Use parametric analyses if neither test rejected, otherwise non-parametric
analyses.
 Shapiro-Wilks Shapiro-Wilks
                              Levenes
                                         Levenes Conclusion
                              Test Stat P-value
   Test Stat P-value
     0 977
                 0.890
                               0.425
                                         0.788 USE PARAMETRIC TESTS
*******************
BASIC SUMMARY STATISTICS
Level N
            Mean StdDev
                              StdErr Coef of Var
                                                      95% Conf.Interval
                              0.27 2.90
0.24 2.72
            18.53 0.54
                                                      17.67, 19.38
 Ctrl 4
           18.00 0.49 0.24
18.33 0.56 0.28
18.13 0.21 0.10
18.40 0.29 0.15
 Dose1 4 18.00
                                                        17.22,
                                                                 18.78
                                                       17.44, 19.21
 Dose2 4
                                          3.03
 Dose3 4
Dose4 4
                                          1.14
                                                       17.80,
                                                                18.45
                                          1.60
                                                       17.93,
                                                                18.87
                   Min Max %of Control(means) %Reduction(means)
17.80 19.10 . .
17.40 18.60 97.17 2.83
        Median
Level
 Ctrl
            18.60
            18.00
 Dose1
                                         98.92
 Dose2
            18.20
                    17.80
                              19.10
                                                            1.08
            18.15
                    17.90
                               18.30
                                         97.84
                                                            2.16
 Dose3
 Dose4
            18.40
                    18.10
                              18.70
                                          99.33
                                                             0.67
*******************
PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
   Analysis of Variance (ANOVA) - overall F-test
    Numerator df Denominator df F-stat
                                                P-value
                    15
                                    0.92
                                                0.476
Dunnett - testing each trt mean signif. different than control
Williams - test assumes dose-response relationship, testing negative trend
Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC
      Mean Dunnett Isotonic Williams
                                                       Tukey p-values
Level
                               p-value Dose1 Dose2 Dose3 Dose4
              p-value
                      mean
                                                                        Dose5
       18.53
Ctrl
                       18.53
      18.00 0.305 18.21 0.198
Dose1
                                                 .
Dose2 18.33 0.915 18.21 0.211 0.831 .
Dose3 18.13 0.528 18.21 0.218 0.994 0.965
Dose4 18.40 0.983 18.21 0.223 0.703 0.999
```

Page 45 of 61

0.223 0.703 0.999 0.898

```
******************
NON-PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
   Kruskal-Wallis test - equality among treatment groups
    Degrees of Freedom TestStat P-value
                          3.39
                                    0.494
MannWhit - testing each trt median signif. different from control
Jonckheere - test assumes dose-response relationship, testing negative trend
Level
        Median
                        MannWhit p-value
                                                    Jonckheere p-value
         18.60
 Ctrl
          18.00
                                   0.332
                                                          0.117
  Dose1
  Dose2
           18.20
                                  0.673
                                                          0.328
          18.15
                                   0.341
                                                          0.254
  Dose3
          18.40
                                   0.887
                                                          0.580
  Dose4
 DECREASING TREND TEST SUMMARY LOWEST CONCENTRATION SIGNIF. LESS THAN CONTROL
  Williams
                                           >highest dose (no sign. differences)
                                            >highest dose (no sign. differences)
  Jonckheere
************************
PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
   Analysis of Variance (ANOVA) - overall F-test
    Numerator df Denominator df F-stat
                                                  P-value
                                     0.92
                      15
Dunnett - testing each trt mean signif. different than control
Williams - test assumes dose-response relationship, testing INCREASING trend
Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC
       Mean Dunnett Isotonic Williams
Level
                                                          Tukey p-values
               p-value mean p-value Dose1 Dose2 Dose3 Dose4 Dose5

      Ctrl
      -18.53
      .
      -18.24
      .
      .
      .

      Dosel
      -18.00
      0.305
      -18.24
      0.881
      .
      .

      Dose2
      -18.33
      0.915
      -18.24
      0.905
      0.831
      .

      Dose3
      -18.13
      0.528
      -18.24
      0.917
      0.994
      0.965

Dose4 -18.40 0.983
                         -18.40
                                  0.803
                                            0.703 0.999 0.898
***********************
NON-PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
   3.39
                                     0.494
MannWhit - testing each trt median signif. different from control
Jonckheere - test assumes dose-response relationship, testing INCREASING trend
                        MannWhit p-value
Level
        Median
                                                   Jonckheere p-value
 Ctrl
         -18.60
                                                         0.883
 Dose1
         -18.00
                                  0.332
                                                         0.672
 Dose2 -18.20
                                  0.673
                                                          0.746
 Dose3
          -18.15
                                   0.341
```

Page 46 of 61

Dose4	-18	. 40		0.887			0.420		
INCREAS	SING TRE	END TEST S	UMMARY	LOWEST	CONCENTRA	rion sig	NIF. GREA	TER THAI	N
CONTROL Willi Jonck	ams theere						(no sign.		
			h screen s IABLE VAR(			length	(mm))		
			PARAMETRI	,		- 3	, , ,		
Shapiro- Levenes	Wilks t test fo	test for N or homogen	ormality of vaif neither	of Residua ariance(ab	ils alph solute res	siduals)	alpha		0.05
Shapir	o-Wilks		-Wilks		Levenes		usion		
	Stat 949	P-val 0.35		Test Stat 4.063	P-value 0.020		ON-PARAME	TRIC TE	STS
*****	****	*****	******	*****	*****	*****	*****	*****	* * * *
BASIC SU Level		STATISTICS Mean St	dDev	StdErr	Coef of V	Jar 0	5% Conf.I	nterval	
			0.19	0.10	8.15	var 9	2.05,		
Dosel			0.05	0.03	2.30		2.10,		
Dose2			0.26	0.13	11.76		1 02	2.67	
Dosez	4		0.20	0.13	3.71		1.83, 2.07,	2.33	
Dose3									
DOSE4	4	2.23	0.05	0.02	2.25		2.15,	2.30	
Level	Me	edian	Min	Max %c	of Control	(means)	%Reduct	ion(mean	ns)
Ctrl			2.10	2.50	•		•		
Dose1		2.20	2.10	2.20	92.55		7.4	5	
Dose2		2.20	2.00	2.60	95.74		4.2	26	
Dose3		2.20	2.10	2.30	93.62		6.3	8	
Dose4		2.20	2.20	2.30	94.68		5.3	32	
*****	*****	*****	*******	*****	******	*****	******	*****	* * * *
			use alpha			l tests			
			(ANOVA) -						
Num	erator		ominator d			P-value			
	4	1	5	0.7	7 (	0.561			
Williams	- test	c assumes	rt mean si dose-respo	nse relat	ionship, t	testing :	negative		
Tukey -	two-sic	ded tests,	all possi	lble compa	risons, no	ot used	for NOEC	or LOEC	
Level	Mean	Dunnett	Isotonic	Williams	}		Tukey p-v	ralues	
		p-value	mean	p-value	e Dosel	Dose2	Dose3	Dose4	Dose5
Ctrl	2.35		2.35						
Dose1	2.18	0.343	2.21	0.134					
Dose2	2.25		2.21	0.143	0.956				
Dose3	2.20	0.471	2.21	0.147	0.999	0.990			
Dose4	2.23	0.618	2.21	0.150	0.990	0.999	0.999		
	3			0				•	·

Page 47 of 61

```
*************************
                     - use alpha-level=0.05 for all tests
NON-PARAMETRIC ANALYSES
   Kruskal-Wallis test - equality among treatment groups
    Degrees of Freedom TestStat P-value
                        2.48
MannWhit - testing each trt median signif. different from control
Jonckheere - test assumes dose-response relationship, testing negative trend
                     MannWhit p-value
Level
       Median
                                             Jonckheere p-value
 Ctrl
        2.40
                              0.270
 Dose1
          2.20
                                                  0.089
         2.20
                              0.673
                                                  0.185
 Dose2
          2.20
                              0.335
                                                  0.207
 Dose3
          2.20
                              0.399
                                                  0.326
 Dose4
DECREASING TREND TEST SUMMARY
                            LOWEST CONCENTRATION SIGNIF. LESS THAN CONTROL
  Williams
                                      >highest dose (no sign. differences)
  Jonckheere
                                      >highest dose (no sign. differences)
*************************
PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
   Analysis of Variance (ANOVA) - overall F-test
    Numerator df Denominator df F-stat
                                           P-value
         4
                   15
                                 0.77
                                           0.561
Dunnett - testing each trt mean signif. different than control
Williams - test assumes dose-response relationship, testing INCREASING trend
Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC
Level Mean Dunnett Isotonic Williams
                                                  Tukey p-values
            p-value mean p-value Dose1 Dose2 Dose3 Dose4
                                                                Dose5
                    Ctrl
      -2.35
Dose1 -2.18 0.343
Dose2 -2.25 0.768
       -2.20 0.471
                                            0.990
Dose 3
      -2.23 0.618
                      -2.24
                              0.939
                                      0.990 0.999 0.999
Dose4
************************
NON-PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
   Kruskal-Wallis test - equality among treatment groups
    Degrees of Freedom TestStat P-value
                       2.48
                                0.648
MannWhit - testing each trt median signif. different from control
Jonckheere - test assumes dose-response relationship, testing INCREASING trend
Level
       Median
                     MannWhit p-value
                                            Jonckheere p-value
 Ctrl
        -2.40
         -2.20
                              0.270
                                                  0.911
 Dose1
         -2.20
                              0.673
                                                  0.815
 Dose2
```

Page 48 of 61

Dose3	-2.	20		0.335		0	.793		
Dose4	-2.			0.399			0.674		
INCREAS CONTROL	ING TRE	ND TEST	SUMMARY	LOWEST	CONCENTRAT	ION SIGN	IIF. GREA	TER THAN	4
Willi	ams				>highes	t dose (	no sign.	differe	ences)
Jonck	heere						no sign.		
	1 '1		1		4 7				
	_		rph screen ARIABLE VA		4-d norm hind-	limb )			
			OR PARAMETI	-					
					als alpha				
					solute res				J.05
use para analyses		analyse	s ii neitne	er test rej	jected, othe	erwise n	ion-param	etric	
-		Shani	ro-Wilks	Levenes	Levenes	Conclu	igi on		
_	Stat	_	alue	Test Stat		COIICI	151011		
	, , , , , , , , , , , , , , , , , , , ,		4140			NO DAT	A FOR TE	ST	
*****	*****	*****	*****	*****	*****	*****	******	*****	***
BASIC SU	MMARY S	TATISTI	CS						
Level			StdDev	StdErr	Coef of Va	ar 95	% Conf.I	nterval	
Ctrl			0.00	0.00	0.00		. ,	•	
Dose1		0.10		0.00	0.00		. ,	•	
Dose2	4	0.10	0.00	0.00	0.00		• ,	•	
Dose3			0.00	0.00	0.00		. ,	•	
Dose4	4	0.10	0.00	0.00	0.00		• ,	•	
Level	Me	dian	Min	Max %c	of Control(	meang)	&Peduat	ion(mear	ag )
Ctrl	1410	0.10	0.10	0.10	or concrot(	.licaris /	*Reduct	TOII ( IIIE a I	.15 /
Dosel		0.10	0.10	0.10	100.00		0.0	0	
Dose2		0.10	0.10	0.10	100.00		0.0		
Dose3		0.10	0.10	0.10	100.00		0.0		
Dose4		0.10	0.10	0.10	100.00		0.0	0	
*****	*****	*****	*****	*****	******	*****	*****	*****	****
PARAMETR					05 for all	tests			
	-		ce (ANOVA)						
Nun	merator	df D	enominator	df F-st	at P	-value			
•	•	•	1 .	•	•	•	•	٠	•
Dunnott	toati	ng oagh	trt moon	nianif dif	ferent than	n gontro	. 7		
		_		_	ionship, te			trand	
					risons, no				
Tukey -	CWO-SIO	ieu test	s, all poss	sible compa	ilisons, no	t useu I	OI NOEC	OI HOEC	
Level	Mean	Dunnet	t Isotonio	c Williams		т	ukey p-v	alues	
ПСАСТ	rican	p-valu		p-value		Dose2	Dose3	Dose4	Dose5
		r vara		F varac	20001		20203	20201	20003
Ctrl	0.10								
Dose1	0.10		•						
Dose2	0.10		•	•					

Page 49 of 61

Dose3	0.10								
Dose4	0.10	•		•					
******	*****	*****	*****					*****	****
		ANALYSES		lpha-level			sts		
			- equality	_	_	roups			
Deg		Freedom		P-valu					
	4		0.00	1.00	U				
MannWhit	- test	ing each	trt median	gianif d	ifferent	from co	ntrol		
			s dose-res					ve trend	
0 0110111100			2 4020 102	Police role	010110111	, 0000111	J 110Jaor		
Level	Medi	an	MannWhi	t p-value		Jonckh	eere p-va	alue	
Ctrl		10							
Dosel	0.	10		1.000					
Dose2				1.000					
Dose3				1.000			•		
Dose4	0.	10		1.000			•		
DECDEAC	ITMO MDE	IND BEIOR O	TIMBAD D37	TOMEGE C		mton ata	NIE IEG		ONTERDOT
DECREAS Willi	_	ND IEST S	UMMARY	LOWEST C	Dose1	IIION SIG	NIF. LES	S THAN C	ONIROL
	heere				Dosel Dosel				
0 OIICh	IIICCIC				DODCI				
*****	*****	*****	*****	*****	*****	*****	*****	*****	***
PARAMETR	RIC ANAL	YSES -	use alpha	-level=0.0	5 for al	l tests			
Anal	ysis of		(ANOVA) -						
Num	nerator	df Den	ominator d	f F-sta	.t	P-value			
•	•	. 1	•	•	•	•	•	•	•
Dunnott	toati	ng oagh t	rt mean si	anif diff	oront th	an contr	o 1		
		_	dose-respo	_				NG trend	
			all possi						
2		,	1		,				
Level	Mean	Dunnett	Isotonic	Williams			Tukey p-	values	
		p-value	mean	p-value	Dose1	Dose2	Dose3	Dose4	Dose5
Ctrl	-0.10	•	•	•	•	•	•	•	•
	-0.10	•	•	•	•	•	•	•	•
Dose2	-0.10	•	•	•	•	•	•	•	•
Dose3 Dose4	-0.10 -0.10	•	•	•	•	•	•	•	•
DOSE4	-0.10	•	•	•	•	•	•	•	•
*****	*****	*****	*****	*****	*****	*****	*****	*****	***
NON-PARA	METRIC	ANALYSES	- use a	lpha-level	=0.05 fo	r all te	sts		
			- equality	-					
		Freedom	TestStat			-			
	4		0.00	1.00	0				
			trt median						
Jonckhee	ere - te	st assume	s dose-res	ponse rela	tionship	, testin	g INCREA	SING tre	nd
								,	
Level	Medi	an	MannWhi	t p-value		Jonckh	eere p-va	aıue	

Page 50 of 61

							LY MUUD	Nullibel 40	3317002
Ctrl	-0	.10							
Dose1	-0	.10		1.000					
Dose2		.10		1.000					
	-0	.10		1.000					
	-0			1.000					
INCREA CONTROL		END TEST	SUMMARY	LOWEST	CONCENTRA	TION SIG	NIF. GREA	ATER THAN	1
Will	iams				Dose1				
Jonc	kheere				Dose1				
L		1			4 -1				
	_		rph screen ARIABLE VAR	_		edian))			
			OR PARAMETR						
			Normality						
			eneity of v						0.05
_		analyses	s if neithe	r test rej	ected, ot.	herwise :	non-para	metric	
analyse		Gl		T	T	G 7			
			ro-Wilks				usion		
	.944		alue 284	6.750			OM_DADAMI	בייסדר ייניס	TTC
U	. 244	0.2	204	0.750	0.003	USE IV	ON-PARAM	EIRIC IES	010
*****	*****	*****	*****	*****	*****	*****	*****	*****	***
BASIC S	UMMARY	STATISTI	CS.						
			StdDev	StdErr	Coef of '	Var 9	5% Conf.	Interval	
		58.50		0.29	0.99		57.58,		
					0.85		57.95,		
Dose2	4	58.75 59.00	0.00	0.25 0.00	0.00		. ,		
Dose3	4	58.50	0.58	0.29	0.99		. , 57.58,	59.42	
	4			0.25			57.45,	59.05	
		edian			f Control	(means)	%Reduct	tion(mean	ıs)
Ctrl		58.50		59.00	•		•		
Dose1		59.00	58.00	59.00	100.43		-0.4		
Dose2			59.00	59.00			-0.8		
Dose3		58.50		59.00	100.00		0.0		
Dose4		58.00	58.00	59.00	99.57		0.4	43	
*****	******	*****	*****	******	******	*****	*****	*******	***
			- use alph						
			ce (ANOVA)						
	merator		enominator			P-value			
Nu	Merator 4	ar De	15	1.3		0.284			
	7		1.0	1.3		0.201			
Dunnett	- test	ing each	trt mean s	ianif. dif	ferent th	an contr	ol		
			s dose-resp					trend	
			s, all poss						
	23 21		, FODD	,	, 11				
Level	Mean	Dunnett	Isotonic	Williams			Tukey p-	values	
		p-value		p-value		Dose2	Dose3	Dose4	Dose5
Ctrl	58.50		58.75					•	

Page 51 of 61

```
58.75
58.75
                             0.841
0.869
      58.75 0.874
Dose1
                                     0.946
      59.00 0.421
Dose2
Dose3
       58.50 1.000
                     58.50
                               0.636
                                       0.946
                                              0.599
      58.25 0.874
                      58.25
                               0.322
                                       0.599
                                             0.234
                                                     0.946
Dose4
********************
NON-PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
   Kruskal-Wallis test - equality among treatment groups
    Degrees of Freedom TestStat P-value
                        5.15
                                 0.273
MannWhit - testing each trt median signif. different from control
Jonckheere - test assumes dose-response relationship, testing negative trend
Level
        Median
                      MannWhit p-value
                                              Jonckheere p-value
         58.50
 Ctrl
                                                    0.753
 Dose1
         59.00
                               0.624
                                                    0.941
         59.00
                               0.223
 Dose2
                              1.000
                                                    0.592
 Dose3
         58.50
         58.00
                               0.624
 Dose4
                                                    0.173
DECREASING TREND TEST SUMMARY LOWEST CONCENTRATION SIGNIF. LESS THAN CONTROL
  Williams
                                       >highest dose (no sign. differences)
  Jonckheere
                                       >highest dose (no sign. differences)
************************
PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
   Analysis of Variance (ANOVA) - overall F-test
    Numerator df Denominator df F-stat
                                            P-value
                   15
        4
                                  1.39
                                            0.284
Dunnett - testing each trt mean signif. different than control
Williams - test assumes dose-response relationship, testing INCREASING trend
Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC
       Mean Dunnett Isotonic Williams
Level
                                                   Tukey p-values
            p-value mean p-value Dose1 Dose2 Dose3 Dose4
                                                                   Dose5
                      -58.50
      -58.50
Ctrl
                            0.427
             0.874
Dose1
      -58.75
                      -58.63
      -59.00 0.421
                      -58.63
                                       0.946
Dose2
                               0.455
Dose3
      -58.50
              1.000
                      -58.63
                                0.471
                                       0.946
                                              0.599
      -58.25 0.874
                                             0.234
Dose4
                      -58.63
                               0.481
                                       0.599
                                                     0.946
*********************
NON-PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
   Kruskal-Wallis test - equality among treatment groups
    Degrees of Freedom TestStat P-value
                        5.15
                                 0.273
MannWhit - testing each trt median signif. different from control
```

Page 52 of 61

Jonckheere - test assumes dose-response relationship, testing INCREASING trend

Williams										
Ctrl	Level	Medi	an	MannWh:	it p-va	lue	Jonckl	neere p-va	lue	
Dose2	Ctrl	-58.	50							
Dose2	Dose1	-59.	00		0.62	4		0.247		
Dose3	Dose2				0.22	3		0.059		
Dose4								0.408		
Williams										
Williams	INCREAS	ING TRE	ND TEST	SUMMARY	LOWE	ST CONCENTRAT:	ION SIG	GNIF. GREA	TER THAN	1
test for amphib metamorph screen study - 2,4-d ANALYSIS RESULTS FOR VARIABLE VAR07 ( 21-d wet weight (g) )  TESTS OF ASSUMPTIONS FOR PARAMETRIC ANALYSIS Shapiro-Wilks test for Normality of Residuals alpha-level=0.01 Levenes test for homogeneity of variance(absolute residuals) alpha-level=0.05 Use parametric analyses if neither test rejected, otherwise non-parametric analyses. Shapiro-Wilks Shapiro-Wilks Levenes Levenes Conclusion Test Stat P-value Test Stat P-value 0.948 0.345 1.312 0.310 USE PARAMETRIC TESTS  *********************************	CONTROL							,	11.55	,
Test for amphib metamorph screen study - 2,4-d ANALYSIS RESULTS FOR VARIABLE VARO7 ( 21-d wet weight (g) )  TESTS OF ASSUMPTIONS FOR PARAMETRIC ANALYSIS Shapiro-Wilks test for Normality of Residuals alpha-level=0.01 Levenes test for homogeneity of variance(absolute residuals) alpha-level=0.05 Use parametric analyses if neither test rejected, otherwise non-parametric analyses.  Shapiro-Wilks Shapiro-Wilks Levenes Levenes Conclusion Test Stat P-value Test Stat P-value 0.948 0.345 1.312 0.310 USE PARAMETRIC TESTS  *********************************										
ANALYSIS RESULTS FOR VARIABLE VARO7 ( 21-d wet weight (g) )  TESTS OF ASSUMPTIONS FOR PARAMETRIC ANALYSIS Shapiro-Wilks test for homogeneity of variance(absolute residuals) alpha-level=0.05 Levenes test for homogeneity of variance(absolute residuals) alpha-level=0.05 Use parametric analyses if neither test rejected, otherwise non-parametric analyses.  Shapiro-Wilks Shapiro-Wilks Levenes Levenes Conclusion Test Stat P-value Test Stat P-value 0.948 0.345 1.312 0.310 USE PARAMETRIC TESTS  *********************************	Joncki	heere				>highes	t dose	(no sign.	differe	ences)
Shapiro-Wilks test for Normality of Residuals alpha-level=0.01		_		-	-	·	t (g)	)		
Shapiro-Wilks test for Normality of Residuals alpha-level=0.01	TECTC ∩E	7 CCIIMD	TTONC E	OD DADAMETED:	T ANTAT	VCTC				
Levenes test for homogeneity of variance(absolute residuals) alpha-level=0.05							a-leve	1=0 01		
Use parametric analyses if neither test rejected, otherwise non-parametric analyses.  Shapiro-Wilks Shapiro-Wilks Levenes Levenes Conclusion  Test Stat P-value Test Stat P-value 0.948 0.345 1.312 0.310 USE PARAMETRIC TESTS  *********************************									-level=f	0.05
Analyses.  Shapiro-Wilks Shapiro-Wilks Test Stat P-value 0.948 0.345 1.312 0.310 USE PARAMETRIC TESTS  *********************************										
Shapiro-Wilks   Shapiro-Wilks   Test Stat   P-value   Test Stat   Test   Test Stat   P-value   Test Stat   Test   Test Stat   P-value   Test Stat   Test   Test	_									
Test Stat			Shapi	ro-Wilks	Levene	s Levenes	Conc	lusion		
### BASIC SUMMARY STATISTICS  Level N Mean StdDev StdErr Coef of Var 95% Conf.Interval Ctrl 4 1.63 0.17 0.08 10.32 1.36, 1.89   ### Dosel 4 1.73 0.12 0.06 6.75 1.54, 1.91   ### Dose2 4 1.61 0.24 0.12 14.76 1.23, 1.99   ### Dose3 4 1.75 0.06 0.03 3.23 1.66, 1.84   ### Dose4 4 1.91 0.22 0.11 11.42 1.56, 2.26    Level Median Min Max % Control(means) % Reduction(means) Ctrl 1.64 1.43 1.78			P-v	alue :						
BASIC SUMMARY STATISTICS   Level N   Mean   StdDev   StdErr   Coef of Var   95% Conf.Interval   Ctrl   4   1.63   0.17   0.08   10.32   1.36,   1.89   Dosel   4   1.73   0.12   0.06   6.75   1.54,   1.91   Dose2   4   1.61   0.24   0.12   14.76   1.23,   1.99   Dose3   4   1.75   0.06   0.03   3.23   1.66,   1.84   Dose4   4   1.91   0.22   0.11   11.42   1.56,   2.26							USE I	PARAMETRIC	TESTS	
BASIC SUMMARY STATISTICS   Level N   Mean   StdDev   StdErr   Coef of Var   95% Conf.Interval   Ctrl   4   1.63   0.17   0.08   10.32   1.36   1.89   Dosel   4   1.73   0.12   0.06   6.75   1.54   1.91   Dose2   4   1.61   0.24   0.12   14.76   1.23   1.99   Dose3   4   1.75   0.06   0.03   3.23   1.66   1.84   Dose4   4   1.91   0.22   0.11   11.42   1.56   2.26										
Level N Mean StdDev StdErr Coef of Var 95% Conf.Interval Ctrl 4 1.63 0.17 0.08 10.32 1.36, 1.89  Dosel 4 1.73 0.12 0.06 6.75 1.54, 1.91  Dose2 4 1.61 0.24 0.12 14.76 1.23, 1.99  Dose3 4 1.75 0.06 0.03 3.23 1.66, 1.84  Dose4 4 1.91 0.22 0.11 11.42 1.56, 2.26  Level Median Min Max % of Control(means) % Reduction(means)  Ctrl 1.64 1.43 1.78  Dosel 1.76 1.57 1.82 106.27 -6.27  Dose2 1.69 1.26 1.80 99.17 0.83  Dose3 1.75 1.70 1.82 107.92 -7.92  Dose4 1.93 1.63 2.15 117.51 -17.51  ***********************************	*****	*****	*****	*****	*****	*****	*****	*****	******	***
Ctrl 4 1.63 0.17 0.08 10.32 1.36, 1.89  Dosel 4 1.73 0.12 0.06 6.75 1.54, 1.91  Dose2 4 1.61 0.24 0.12 14.76 1.23, 1.99  Dose3 4 1.75 0.06 0.03 3.23 1.66, 1.84  Dose4 4 1.91 0.22 0.11 11.42 1.56, 2.26  Level Median Min Max % of Control(means) %Reduction(means)  Ctrl 1.64 1.43 1.78  Dose1 1.76 1.57 1.82 106.27 -6.27  Dose2 1.69 1.26 1.80 99.17 0.83  Dose3 1.75 1.70 1.82 107.92 -7.92  Dose4 1.93 1.63 2.15 117.51 -17.51  ***********************************	BASIC SUN	MMARY S	TATISTI	CS						
Dosel 4 1.73 0.12 0.06 6.75 1.54, 1.91 Dose2 4 1.61 0.24 0.12 14.76 1.23, 1.99 Dose3 4 1.75 0.06 0.03 3.23 1.66, 1.84 Dose4 4 1.91 0.22 0.11 11.42 1.56, 2.26  Level Median Min Max % Control(means) % Reduction(means) Ctrl 1.64 1.43 1.78 Dose1 1.76 1.57 1.82 106.27 -6.27 Dose2 1.69 1.26 1.80 99.17 0.83 Dose3 1.75 1.70 1.82 107.92 -7.92 Dose4 1.93 1.63 2.15 117.51 -17.51  ***********************************					StdErr	Coef of Va	ar 🤉			
Dose2 4 1.61 0.24 0.12 14.76 1.23, 1.99 Dose3 4 1.75 0.06 0.03 3.23 1.66, 1.84 Dose4 4 1.91 0.22 0.11 11.42 1.56, 2.26  Level Median Min Max % of Control(means) % Reduction(means) Ctrl 1.64 1.43 1.78 . Dose1 1.76 1.57 1.82 106.27 -6.27 Dose2 1.69 1.26 1.80 99.17 0.83 Dose3 1.75 1.70 1.82 107.92 -7.92 Dose4 1.93 1.63 2.15 117.51 -17.51  ***********************************	Ctrl	4						1.36,	1.89	
Dose3 4 1.75 0.06 0.03 3.23 1.66, 1.84 Dose4 4 1.91 0.22 0.11 11.42 1.56, 2.26  Level Median Min Max % Control(means) % Reduction(means) Ctrl 1.64 1.43 1.78 Dose1 1.76 1.57 1.82 106.27 -6.27 Dose2 1.69 1.26 1.80 99.17 0.83 Dose3 1.75 1.70 1.82 107.92 -7.92 Dose4 1.93 1.63 2.15 117.51 -17.51  ***********************************		4	1.73	0.12				1.54,	1.91	
Dose4 4 1.75 0.06 0.03 3.23 1.66, 1.84  Dose4 4 1.91 0.22 0.11 11.42 1.56, 2.26  Level Median Min Max % Control(means) % Reduction(means)  Ctrl 1.64 1.43 1.78  Dose1 1.76 1.57 1.82 106.27 -6.27  Dose2 1.69 1.26 1.80 99.17 0.83  Dose3 1.75 1.70 1.82 107.92 -7.92  Dose4 1.93 1.63 2.15 117.51 -17.51  ***********************************		4						1.23,	1.99	
Level Median Min Max % of Control(means) % Reduction(means) Ctrl 1.64 1.43 1.78  Dosel 1.76 1.57 1.82 106.27 -6.27  Dose2 1.69 1.26 1.80 99.17 0.83  Dose3 1.75 1.70 1.82 107.92 -7.92  Dose4 1.93 1.63 2.15 117.51 -17.51  ***********************************		4	1./5							
Ctrl 1.64 1.43 1.78	Dose4	4	1.91	0.22	0.11	11.42		1.56,	2.26	
Dosel 1.76 1.57 1.82 106.27 -6.27 Dose2 1.69 1.26 1.80 99.17 0.83 Dose3 1.75 1.70 1.82 107.92 -7.92 Dose4 1.93 1.63 2.15 117.51 -17.51  ***********************************	Level	Me	dian	Min	Max	%of Control(	means)	%Reduct	ion(mear	ns)
Dose2 1.69 1.26 1.80 99.17 0.83 Dose3 1.75 1.70 1.82 107.92 -7.92 Dose4 1.93 1.63 2.15 117.51 -17.51  ***********************************	Ctrl		1.64	1.43	1.78					
Dose3 1.75 1.70 1.82 107.92 -7.92 Dose4 1.93 1.63 2.15 117.51 -17.51  ***********************************	Dose1				1.82	106.27		-6.2	7	
Dose4 1.93 1.63 2.15 117.51 -17.51  ***********************************	Dose2		1.69	1.26	1.80	99.17		0.8	3	
**************************************	Dose3		1.75	1.70	1.82	107.92		-7.9	2	
PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests Analysis of Variance (ANOVA) - overall F-test Numerator df Denominator df F-stat P-value 4 15 1.94 0.156  Dunnett - testing each trt mean signif. different than control Williams - test assumes dose-response relationship, testing negative trend Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC  Level Mean Dunnett Isotonic Williams Tukey p-values	Dose4		1.93	1.63	2.15	117.51		-17.5	1	
PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests Analysis of Variance (ANOVA) - overall F-test Numerator df Denominator df F-stat P-value 4 15 1.94 0.156  Dunnett - testing each trt mean signif. different than control Williams - test assumes dose-response relationship, testing negative trend Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC  Level Mean Dunnett Isotonic Williams Tukey p-values	*****	*****	******	*****	******	*****	*****	******	******	***
Analysis of Variance (ANOVA) - overall F-test  Numerator df Denominator df F-stat P-value  4 15 1.94 0.156  Dunnett - testing each trt mean signif. different than control  Williams - test assumes dose-response relationship, testing negative trend  Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC  Level Mean Dunnett Isotonic Williams Tukey p-values										
Numerator df Denominator df F-stat P-value 4 15 1.94 0.156  Dunnett - testing each trt mean signif. different than control Williams - test assumes dose-response relationship, testing negative trend Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC  Level Mean Dunnett Isotonic Williams Tukey p-values										
4 15 1.94 0.156  Dunnett - testing each trt mean signif. different than control Williams - test assumes dose-response relationship, testing negative trend Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC  Level Mean Dunnett Isotonic Williams Tukey p-values										
Dunnett - testing each trt mean signif. different than control Williams - test assumes dose-response relationship, testing negative trend Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC Level Mean Dunnett Isotonic Williams Tukey p-values	Name		ar D							
Williams - test assumes dose-response relationship, testing negative trend Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC Level Mean Dunnett Isotonic Williams Tukey p-values		7		13		1.94	.130			
Williams - test assumes dose-response relationship, testing negative trend Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC Level Mean Dunnett Isotonic Williams Tukey p-values	Dunnett -	- testi	ng each	trt mean s	ianif (	different that	n conti	col		
Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC  Level Mean Dunnett Isotonic Williams Tukey p-values									trend	
• •				_		_	_			
* *	-									
p-value mean p-value Dose1 Dose2 Dose3 Dose4 Dose5	Level	Mean	Dunnet	t Isotonic	Willi	ams		Tukey p-v	alues	
			p-valu	e mean	p-va	lue Dosel	Dose2	Dose3	Dose4	Dose5

Page 53 of 61

```
    Ctrl
    1.63
    .
    1.73
    .
    .
    .

    Dosel
    1.73
    0.863
    .
    .
    .

    Dose2
    1.61
    1.000
    1.73
    0.889
    0.875
    .

    Dose3
    1.75
    0.681
    1.73
    0.902
    0.999
    0.770

    Dose4
    1.81
    0.104
    1.73
    0.602
    0.202
    0.202

       1.91 0.104
                                   0.909
Dose4
                          1.73
                                              0.580
                                                    0.157 0.708
********************
NON-PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
   Kruskal-Wallis test - equality among treatment groups
    Degrees of Freedom TestStat P-value
                            5.41
                                      0.247
MannWhit - testing each trt median signif. different from control
Jonckheere - test assumes dose-response relationship, testing negative trend
                          MannWhit p-value
Level
          Median
                                                      Jonckheere p-value
 Ctrl
           1 64
                                    0.346
                                                            0.876
 Dose1
            1.76
                                   1.000
                                                            0.500
 Dose2
           1.69
           1.75
                                    0.494
                                                            0.772
 Dose3
 Dose4
           1.93
                                    0.156
                                                            0.969
DECREASING TREND TEST SUMMARY LOWEST CONCENTRATION SIGNIF. LESS THAN CONTROL
  Williams
                                             >highest dose (no sign. differences)
                                              >highest dose (no sign. differences)
  Jonckheere
******************
                    - use alpha-level=0.05 for all tests
PARAMETRIC ANALYSES
    Analysis of Variance (ANOVA) - overall F-test
    Numerator df Denominator df F-stat
                                                    P-value
                      15
                                                    0.156
                                       1.94
Dunnett - testing each trt mean signif. different than control
Williams - test assumes dose-response relationship, testing INCREASING trend
Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC
        Mean Dunnett Isotonic Williams
Level
                                                            Tukey p-values
               p-value mean p-value Dose1 Dose2 Dose3 Dose4
                                                                             Dose5
                          ...67 0.428
-1.67 0.457
-1.75 0
Ctrl
        -1.63
                       -1.67
-1.67
-1.75
        -1.73 0.819
                                                      .
Dose1
        -1.61 1.000
-1.75 0.681
Dose2
                                              0.875
                                                      0.770
Dose3
                                              0.999
       -1.91 0.104
                                             0.580 0.157
Dose4
                                                             0.708
********************
NON-PARAMETRIC ANALYSES
                         - use alpha-level=0.05 for all tests
   Kruskal-Wallis test - equality among treatment groups
    Degrees of Freedom TestStat P-value
           4
                             5.41
                                       0.247
```

Page 54 of 61

Jonckheere - test assumes dose-response relationship, testing INCREASING trend

MannWhit - testing each trt median signif. different from control

Level	Median	MannWh	nit p-value		Jonckheere p-value	
Ctrl	-1.64					
Dose1	-1.76		0.346		0.124	
Dose2	-1.69		1.000		0.500	
Dose3	-1.75		0.494		0.228	
Dose4	-1.93		0.156		0.031	
INCREAS	ING TREND T	EST SUMMARY	LOWEST	CONCENTRAT	ON SIGNIF. GREATER THAN	
CONTROL						
Willi	ams			Dose4		
Jonck:				Dose4		
test for	amphib met	amorph screen	study - 2.	4-d		
	_	R VARIABLE VAF	_		ength (mm) )	
111111111111111111111111111111111111111	1120210 10	VIII.		211 (0110 1	511 <u>3</u> 611 (	
TESTS OF	ASSUMPTION	S FOR PARAMETE	RTC ANALYST	S		
		for Normality	-	_	a-level=0.01	
					duals) alpha-level=0.	05
					erwise non-parametric	
analyses		ybeb ii neiene	LI CCDC ICJ	ccca, oun	rwibe non parametric	
		apiro-Wilks	T.eveneg	T.evenes	Conclusion	
	Stat		Test Stat	P-value	Coliciusion	
	981	0.945	5.995	0.004	USE NON-PARAMETRIC TEST	'C
0.	701	0.743	3.773	0.004	ODE NON PARAMETRIC TEST	D
*****	*****	*****	*****	*****	*******	**
	MMARY STATI					
Level i			StdErr	Coef of Va	ar 95% Conf.Interval	
Ctrl	4 26.78		0.70	5.21	24.56, 28.99	
Dose1			0.70	2.04	26.51, 28.29	
			0.28			
Dosez	4 26.85 4 27.68	1.03	0.81	6.06 1.12	24.26, 29.44	
Dose3 Dose4	4 28.08		0.15	1.12	27.18, 28.17 27.43, 28.72	
Dose4	4 20.00	0.40	0.20	1.44	27.43, 20.72	
Level	Median	Min	Max %o	f Control(r	means) %Reduction(means	. \
Ctrl	26.80		28.20	r control(i	means) % Reduction ( means	' /
	27.25			102 22		
Dose1 Dose2			28.20 28.70	102.33	-2.33 -0.28	
	26.90			100.28 103.36	-3.36	
Dose3 Dose4	27.60 28.10		28.10 28.50	103.36	-3.36 -4.86	
Dose4	20.10	27.00	20.50	104.00	-4.00	
******	********	*****	********	*******	******	**
						~ ~
		- use alph			tests	
	-	iance (ANOVA)			1	
Num	erator df	Denominator			-value	
	4	15	1.1	8 0	.361	
D		1		F		
		ach trt mean s				
					esting negative trend	
Tukey -	two-sided t	ests, all poss	sible compa	risons, not	used for NOEC or LOEC	
	., -		**** 7 7 1		m 1	
Level		nett Isotonio			Tukey p-values	
	p-v	alue mean	p-value	Dose1	Dose2 Dose3 Dose4	Dose5

Page 55 of 61

```
26.78
Ctrl
                      27.36
Dosel 27.40 0.798 27.36 0.859
Dose2 26.85 1.000 27.36 0.885
Dose3 27.68 0.550 27.36 0.899
                                       0.937
                                             0.779
                                       0.995
      28.08 0.253
                      27.36
                              0.907
                                       0.877 0.460
Dose4
*******************
NON-PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
   Kruskal-Wallis test - equality among treatment groups
    Degrees of Freedom TestStat P-value
                        4.19
                                 0.381
MannWhit - testing each trt median signif. different from control
Jonckheere - test assumes dose-response relationship, testing negative trend
Level
        Median
                      MannWhit p-value
                                              Jonckheere p-value
 Ctrl
         26.80
         27.25
                               0.780
                                                    0.668
 Dose1
         26.90
                              1.000
                                                   0.529
 Dose2
                              0.780
 Dose3
         27.60
                                                    0.787
 Dose4
         28.10
                               0.235
                                                    0.971
DECREASING TREND TEST SUMMARY LOWEST CONCENTRATION SIGNIF. LESS THAN CONTROL
                                        >highest dose (no sign. differences)
  Williams
  Jonckheere
                                        >highest dose (no sign. differences)
******************
PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
   Analysis of Variance (ANOVA) - overall F-test
    Numerator df Denominator df F-stat
                                            P-value
         4
                   15
                                 1.18
                                            0.361
Dunnett - testing each trt mean signif. different than control
Williams - test assumes dose-response relationship, testing INCREASING trend
Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC
Level Mean Dunnett Isotonic Williams
                                                   Tukey p-values
            p-value mean p-value Dose1 Dose2 Dose3 Dose4
                                                                  Dose5
                            0.377
      -26.78
                      -26.78
Ctrl
      -27.40 0.798
Dose1
                      -27.13
             1.000
                      -27.13
                                       0.937
Dose2 -26.85
                               0.403
Dose3 -27.68 0.550
                      -27.68
                               0.150
                                       0.995
                                              0.779
Dose4 -28.08 0.253
                      -28.08
                               0.058
                                       0.877 0.460
                                                     0.979
********************
NON-PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
   Kruskal-Wallis test - equality among treatment groups
    Degrees of Freedom TestStat P-value
         4
                       4.19
                                 0.381
```

Page 56 of 61

MannWhit - testing each trt median signif. different from control

```
Jonckheere - test assumes dose-response relationship, testing INCREASING trend
 Level
          Median
                           MannWhit p-value
                                                      Jonckheere p-value
 Ctrl
          -26.80
          -27.25
                                    0.780
                                                             0.332
  Dose1
  Dose2
          -26.90
                                    1.000
                                                             0.471
                                    0.780
                                                             0.213
  Dose3
          -27.60
                                    0.235
                                                             0.029
 Dose4
          -28.10
 INCREASING TREND TEST SUMMARY
                                   LOWEST CONCENTRATION SIGNIF. GREATER THAN
CONTROL
  Williams
                                               >highest dose (no sign. differences)
   Jonckheere
                                              Dose4
test for amphib metamorph screen study - 2,4-d
ANALYSIS RESULTS FOR VARIABLE VAR09 ( 21-d hind-limb length (mm) )
TESTS OF ASSUMPTIONS FOR PARAMETRIC ANALYSIS
Shapiro-Wilks test for Normality of Residuals -- alpha-level=0.01
Levenes test for homogeneity of variance(absolute residuals) -- alpha-level=0.05
Use parametric analyses if neither test rejected, otherwise non-parametric
analyses.
  Shapiro-Wilks Shapiro-Wilks
                                Levenes
                                            Levenes Conclusion
    Test Stat P-value
                                 Test Stat P-value
     0.955
                   0.443
                                 0.118
                                             0.974 USE PARAMETRIC TESTS
*******************
            StdDev StdErr Coef of Var
17.23 0.77 0.38 4.46
18.00 1.02 0.51
17.08
BASIC SUMMARY STATISTICS
 Level N
                                                           95% Conf.Interval
                                                             16.00, 18.45
 Ctrl
       4
 Dosel 4
                                                             16.38,
                                                                      19.62
             17.08 0.91 0.46
17.48 1.00 0.50
16.20 0.98 0.49
 Dose2 4
                                                             15.62,
                                                                      18.53
 Dose3 4
                                              5.71
                                                             15.89, 19.06
 Dose4 4
                                              6.07
                                                             14.64, 17.76

        Median
        Min
        Max
        % of Control(means)
        % Reduction(means)

        17.20
        16.40
        18.10
        .
        .

        17.70
        17.20
        19.40
        104.50
        -4.50

            Median
 Level
 Ctrl
  Dose1
                      16.20
                                                                  0.87
  Dose2
              16.90
                                  18.30
                                             99.13
              17.60
                       16.30
                                  18.40
                                             101.45
  Dose3
                                                                  -1.45
              16.45
                       14.80
                                  17.10
                                              94.05
                                                                   5.95
*******************
PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
    Analysis of Variance (ANOVA) - overall F-test
     Numerator df Denominator df F-stat
                                                    P-value
                                       1.96
                       15
                                                    0.153
Dunnett - testing each trt mean signif. different than control
Williams - test assumes dose-response relationship, testing negative trend
Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC
Level
       Mean Dunnett Isotonic Williams
                                                             Tukey p-values
```

Page 57 of 61

		1		1	D1	Daga	Daga	D = = = 1	Danaf
		p-value	mean	p-value	Dose1	Dose2	Dose3	Dose4	Dose5
Ctrl	17.23		17.61						
Dose1	18.00	0.607	17.61	0.799		•		•	•
Dose2	17.08		17.28	0.650	0.642	•	•	•	•
Dose3		0.986	17.28	0.669	0.930	0.973	•	•	•
Dose4		0.378	16.20	0.094	0.100	0.686	0.350	•	
*****	*****	****	*****	*****	*****	*****	******	*****	***
NON-PAR	AMETRIC	ANALYSES	- 11SE a	.lpha-level:	=0 05 fo	r all te	sts		
				among trea			565		
		Freedom	TestStat		_	Loaps			
20	4	Trecaom	5.67	0.22	_				
				signif. di ponse relat				wa trand	
UUIICKIIE	ere - ce	sc assumer	duse-les	ponse rela	LIOHSHIP	, cescin	g Hegati	ve crena	
Level Ctrl	Medi 17.		MannWhi	t p-value		Jonckh	eere p-v	alue	
Dose1				0.413			0.845		
Dose1				0.889			0.384		
Dose2				0.678			0.500		
Dose3				0.283			0.086		
20201		10		0.200					
DECREA	SING TRE	ND TEST SU	JMMARY	LOWEST CO	ONCENTRA'	TION SIG	NIF. LES	S THAN C	ONTROL
Will	iams					st dose			
Jone	kheere				>highe	st dose	(no sign	. differ	ences)
*****	*****	****	*****	*****	*****	*****	******	*****	***
PARAMET	RIC ANAL	YSES -	use alpha	-level=0.0!	5 for al	l tests			
			_	overall F					
	merator		ominator d			P-value			
	4	15	5	1.96		0.153			
				gnif. diffe					
				nse relatio					
Тикеу -	two-slo	led tests,	all possi	ble compar	isons, n	ot usea :	for NOEC	or LOEC	
Level	Mean	Dunnett	Isotonic	Williams		ŗ	Tukey p-	values	
		p-value	mean	p-value	Dose1	Dose2	Dose3	Dose4	Dose5
Ctrl	-17.23		-17.20						
Dose1	-18.00	0.607	-17.20	0.602	•				
Dose2	-17.08	0.998	-17.20	0.637	0.642				
	-17.48	0.986	-17.20	0.656	0.930	0.973			
	-16.20	0.378	-17.20	0.668	0.100	0.686	0.350	•	•
*****	****	****	· * * * * * * * * * * * * * * * * * * *	*****	****	****	*****	****	****
		ANALYSES		.lpha-level:					
				among trea			515		
		IIS test - Freedom				roups			
Dе	_	rreedou	TestStat						
	4		5.67	0.22	0				

Page 58 of 61

MannWhit - testing each trt median signif. different from control Jonckheere - test assumes dose-response relationship, testing INCREASING trend Jonckheere p-value Median MannWhit p-value Ctrl -17.20Dose1 -17.700.413 0.155 0.889 0.616 Dose2 -16.90 0.678 0.500 Dose3 -17.60Dose4 -16.450.283 0.914 INCREASING TREND TEST SUMMARY LOWEST CONCENTRATION SIGNIF. GREATER THAN CONTROL Williams >highest dose (no sign. differences) Jonckheere >highest dose (no sign. differences) test for amphib metamorph screen study - 2,4-d ANALYSIS RESULTS FOR VARIABLE VAR10 ( 21-d norm hind-limb ) TESTS OF ASSUMPTIONS FOR PARAMETRIC ANALYSIS Shapiro-Wilks test for Normality of Residuals -- alpha-level=0.01 Levenes test for homogeneity of variance(absolute residuals) -- alpha-level=0.05 Use parametric analyses if neither test rejected, otherwise non-parametric analyses. Shapiro-Wilks Shapiro-Wilks Levenes Levenes Conclusion Test Stat P-value Test Stat P-value < .001 0.750 0.573 USE NON-PARAMETRIC TESTS \* BASIC SUMMARY STATISTICS 
 Level
 N
 Mean
 StdDev
 StdErr
 Coef of Var
 95% Conf.Interval

 Ctrl
 4
 0.68
 0.05
 0.03
 7.41
 0.60, 0.75

 Dosel
 4
 0.65
 0.06
 0.03
 8.88
 0.56, 0.74

 Dosel
 4
 0.65
 0.06
 0.03
 8.88
 0.56, 0.74

 0.65
 0.06

 0.65
 0.06

 0.65
 0.06

 0.58
 0.05

 0.02
 Dose3 4 8.88 0.56, 0.74 Dose4 4 8.70 0.50, 0.65 
 Level
 Median
 Min
 Max
 %of Control(means)
 %Reduction(means)

 Ctrl
 0.70
 0.60
 0.70
 .
 .
 .

 Dosel
 0.65
 0.60
 0.70
 96.30
 3.70

 Dose2
 0.65
 0.60
 0.70
 96.30
 3.70

 Dose3
 0.65
 0.60
 0.70
 96.30
 3.70

 Dose4
 0.60
 0.50
 0.60
 85.19
 14.81
 Dose4 0.60 0.50 0.60 85.19 \* PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests Analysis of Variance (ANOVA) - overall F-test Numerator df Denominator df F-stat P-value 15 1.92 Dunnett - testing each trt mean signif. different than control

Page 59 of 61

Williams - test assumes dose-response relationship, testing negative trend Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC

Level	Mean	Dunnett	Isotonic	Williams			Tukey p-	values	
		p-value	mean	p-value	Dose1	Dose2	Dose3	Dose4	Dose5
~. J	0.60		0.60						
Ctrl	0.68		0.68		•	•		•	•
Dosel	0.65		0.65	0.315		•	•	•	•
Dose2	0.65	0.914	0.65	0.337	1.000	•	•	•	•
Dose3	0.65	0.914	0.65	0.349	1.000	1.000			
Dose4	0.58	0.066	0.58	0.013	0.341	0.341	0.341	•	•
++++++			++++++++	*****	+++++++	+++++++	++++++	+++++++	++++
				alpha-level:					^ ^ ^ ^
		ANALYSES		among tre			SLS		
				P-value	_	roups			
Deg	_	Freedom							
	4		5.87	0.20	9				
MannWhit	t – test	ing each	trt mediar	n signif. d	ifferent	from co	ntrol		
				sponse rela				ve trend	
0011011110	010			Police rela	01101111	, 5555211	J 110Jaor		
Level	Medi	ian	MannWhi	t p-value		Jonckh	eere p-v	alue	
Ctrl	0.	. 70							
Dose1	0.	. 65		0.624			0.247		
Dose2		. 65		0.624			0.246		
Dose3		. 65		0.624			0.256		
Dose4		. 60		0.100			0.019		
DOBCI	0.	. 00		0.100			0.019		
DECREAS	SING TRE	END TEST S	UMMARY	LOWEST CO	ONCENTRAT	TION SIG	NIF. LES	S THAN C	ONTROL
Will:	iams				Dose4				
	kheere				Dose4				
*****	******	*****	*****	*****	*****	*****	*****	*****	****
PARAMETI	RIC ANAI	LYSES -	use alpha	a-level=0.0	5 for all	l tests			
Ana	lysis of	Variance	(ANOVA) -	overall F	-test				
Nur	merator	df Den	ominator d	lf F-sta	t I	P-value			
	4	1	.5	1.92	(	0.160			
Dunnett	- testi	ing each t	rt mean si	gnif. diffe	erent tha	an contr	ol		
Williams	s - test	assumes	dose-respo	nse relation	onship, t	testing	INCREASI	NG trend	
Tukey -	two-sid	ded tests,	all possi	ble compar	isons, no	ot used	for NOEC	or LOEC	
Level	Mean		Isotonic				Tukey p-		
		p-value	mean	p-value	Dose1	Dose2	Dose3	Dose4	Dose5
Ctrl	-0.68		-0.64	•	•	•	•	•	•
Dose1		0.914	-0.64	0.880	•	•	•	•	•
Dose2		0.914	-0.64	0.905	1.000				
Dose3	-0.65	0.914	-0.64	0.916	1.000	1.000		•	
Dose4	-0.58	0.066	-0.64	0.924	0.341	0.341	0.341		
				******				*****	***
		ANALYSES		lpha-level:			STS		
				among trea	_	roups			
Deg	_	Freedom	TestStat						
	4		5.87	0.20	9				

Page 60 of 61

MannWhit - testing each trt median signif. different from control Jonckheere - test assumes dose-response relationship, testing INCREASING trend

_			_	
Level	Median	MannWhit	p-value	e Jonckheere p-value
Ctrl	-0.70			•
Dose1	-0.65		0.624	0.753
Dose2	-0.65		0.624	0.754
Dose3	-0.65		0.624	0.744
Dose4	-0.60		0.100	0.981
INCREASI	NG TREND	TEST SUMMARY	LOWEST	CONCENTRATION SIGNIF. GREATER THAN
CONTROL				
William	ms			>highest dose (no sign. differences)
Jonckheere >highest dose (			>highest dose (no sign. differences)	

# DATA EVALUATION RECORD

2,4-DICHLOROPHENOXYACETIC ACID (2,4-D)

Study Type: OCSPP 890.1150, Androgen Receptor Binding (Rat Prostate Cytosol)

EPA Contract No. EP10H001452 Task Assignment No. 2-26-2012 (MRID 48614301)

Prepared for
Health Effects Division
Office of Pesticide Programs
U.S. Environmental Protection Agency
2777 South Crystal Drive
Arlington, VA 22202

Prepared by CSS-Dynamac Corporation 1910 Sedwick Road, Building 100, Suite B Durham, NC 27713

Primary Reviewer:	Signature:	Sendia Nasting
Sandra Hastings	Date:	3/19/2012
Secondary Reviewer:	Signature:	Lectenberg
Scott D. Studenberg, Ph. D., D.A.B.T.	Date:	3/26/2012
	~	Jack Q. Eusy
Program Manager:	Signature: _	
Jack D. Early, M.S.	Date:	3/26/2012
		Jack D. Ewy
Quality Assurance:	Signature:	
Jack D. Early, M.S.	Date:	3/26/2012

This Data Evaluation Record may have been altered by the Health Effects Division subsequent to signing by CSS-Dynamac Corporation personnel.

The US EPA Endocrine Disruptor Screening Program (EDSP) Tier 1 screening battery is comprised of eleven screening assays intended to identify a chemical's likely endocrine bioactivity, i.e., its potential to interact with the estrogen, androgen, or thyroid (E, A, or T) pathways. The robustness of the Tier 1 battery is based on the strengths of each individual assay to identify potential endocrine bioactivity with complementary endpoints within the assay, where available, and redundancy across the battery. Thus, the results of each individual assay should not be considered in isolation but rather should be considered in the context of other assays in the battery as well as Other Scientifically Relevant Information (OSRI). In order to determine if a chemical has the potential to interact with the E, A or T pathways, a Weight of Evidence (WoE) evaluation of Tier 1 assay results, in combination with the findings in the OSRI, should be undertaken (refer to the WoE Document).

Primary Reviewer: Chester Rodriguez, Ph.D.

**Health Effects Division** 

Secondary Reviewer: Greg Akerman, Ph.D.

**Health Effects Division** 

Signature:

Date:

Signature:

Date:

Template version 08/2011

# DATA EVALUATION RECORD

STUDY TYPE: Androgen Receptor Binding (Rat Prostate Cytosol); OCSPP 890.1150

PC CODE: 030001 DP BARCODE: D398638

TXR#: 0052104 CAS No.: 94-75-7

**TEST MATERIAL (PURITY)**: 2,4-D (98.5% a.i.)

**SYNONYMS**: (2,4-dichlorophenoxy) acetic acid

CITATION: LeBaron, M.J., Schisler, M.R., and Visconti, N.R. (2011). Evaluation of 2,4-

Dichlorophenoxy Acetic Acid (2,4-D) In An *In Vitro* Androgen Receptor Binding Assay. The Dow Chemical Company, Toxicology & Environmental Research and Consulting, Midland, MI. Laboratory Project Study ID: 111111, October 27,

2011. MRID #48614301. Unpublished.

SPONSOR: Industry Task Force II on 2,4-D Research Data, c/o McKenna Long & Aldridge

LLP, 1900 K Street NW, Washington, D.C.

**TEST ORDER #**: CON-030001-1

EXECUTIVE SUMMARY: In an androgen receptor (AR) binding assay (MRID 48614301), ventral prostate cytosol from Sprague Dawley rats was used as the source of AR to conduct a competitive binding experiment to measure the binding of a single concentration of [³H]-R1881 (1 nM) in the presence of increasing concentrations (10<sup>-11</sup> to 10<sup>-4</sup> M) of 2,4-D (98.5% purity). The test guideline recommends testing up to 10<sup>-3</sup> M; however, the sponsor selected 10<sup>-4</sup> M as the highest concentration based on *in vivo* toxicokinetic analyses in the rat; concentrations higher than 10<sup>-4</sup> M were not relevant for testing in this assay as they are substantially above the inflection point for linear toxicokinetics (see Appendix B of study report). Ethanol was used as a solvent at a final concentration of <3%. A total of three independent runs were performed, and the assay included dexamethasone as a weak positive control, and R1881 as the ligand reference standard.

Saturation binding experiments were conducted to demonstrate that the AR in the rat prostate cytosol was present at adequate levels and functioning with appropriate affinity for the radiolabeled ligand. The saturation binding experiment resulted in a maximum binding capacity ( $B_{max}$ ) of 3.245 fmol/100 µg protein and the dissociation constant ( $K_d$ ) was 0.4641 nM. Although these values were slightly below the range of values from the validation studies, the results were highly reproducible and all other performance criteria and the competitive binding

assays indicated acceptable performance of the assay. The Scatchard plot indicated a linear response across the concentrations of ligand added. Nonspecific binding as a percent of total binding was less than 20% across the entire concentration range in the saturation binding assays (range 6.2-19.8%, with the exception of the high concentration (10 nM) in one assay, which was 24.6%).

There were no appreciable alterations in R1881 AR binding activity at 2,4-D concentrations ranging from  $10^{-11}$  to  $10^{-4}$  M in the competitive binding experiments, therefore, the log IC<sub>50</sub> and relative binding affinity (RBA) for 2,4-D could not be calculated. The log IC<sub>50</sub> values for R-1881 alone and the positive control, dexamethasone, were -9.0 and -4.4 M, respectively. Compared to R1881, the RBA for dexamethasone was 0.0027 %. In all instances, R1881 and the positive control met the QC performance criteria established in the test guideline.

Based on the results from the three runs, 2,4-D is classified as a Non-binder in the Androgen Receptor Binding Assay.

The study **satisfies** the EDSP Tier 1 Test Order requirements for an Androgen Receptor Binding Assay (OCSPP 890.1150).

**<u>COMPLIANCE</u>**: Signed and dated GLP and Quality Assurance statements were provided in the study report.

# I. MATERIALS AND METHODS

# A. MATERIALS

1. **Test Facility:** The Dow Chemical Company, Toxicology & Environmental

Research and Consulting

**Location:** Midland, MI **Study Director:** Schisler, M.R.

Other Personnel: LeBaron, M.J. (Lead Scientist); Visconti, N.R. (Research Biologist); Gollapudi, B.B.

(Technical reviewer)

**Study Period:** July 11, 2011 – October 27, 2011

2. Test substance: 2,4-D

**Description:** Off-white powder

Source: NuFarm Americas, Inc. (Burr Ridge, IL)

**Lot/Batch #:** 2006 2433 8006-USA (expiry date: March 3, 2013)

**Purity:** 98.5% a.i.

**Solubility:** Up to 315 mg/L in water; soluble in ethanol up to 30 mM

Volatility: 1.9 x 10<sup>-5</sup> Pa at 25°C Stability: 2-yr shelf life Storage conditions: Ambient CAS #: 94-75-7 Molecular weight: 221.0

**Structure:** 

OH

3. Non-labeled ligand: R1881

Supplier: Perkin-Elmer (Boston, MA)

 Catalog #:
 R0908

 Lot #:
 614156

 Purity:
 >97%

 CAS #:
 965-93-5

4. Radioactive ligand: [3H]-R1881

Supplier: Perkin-Elmer (Boston, MA)
Catalog and Batch #: NET590250UC, Lot# 614814

Catalog and Batch #: NET590250UC, Lot# 6148
Date of production: July 1, 2010

**Date of use:** July 11, 2011 to July 28, 2011

Radiochemical purity: >97%
Specific activity: 85.1 Ci/mmol
Concentration of stock: 1.0 mCi/mL

**5. Positive control:** Dexamethasone

Supplier: Sigma (St. Louis, MO)

 Catalog #
 D4902

 Lot #:
 BCBC9269

 Purity:
 98.9%

 CAS #:
 50-02-2

6. Solvent/vehicle control: Ethanol
Justification for choice of None provided

solvent:
Final Concentration: <3%

# B. METHODS

1. <u>Preparation of Rat Ventral Prostate Cytosol</u>: The rat ventral prostate tissue was purchased from Charles River Laboratories (Wilmington, MA). Male Sprague Dawley rats (number not reported) were castrated at approximately 90 days of age and euthanized approximately 24 hours later. The ventral prostate tissues were collected and stored at approximately -80°C until use, and were processed as a batch and used for multiple studies.

The cytosol was prepared by adding low-salt TEDG buffer [0.01 M Tris, 1 mM sodium molybdate, 1.5 mM EDTA, 10% glycerol and 1 mM phenylmethylsulfonyl fluoride (PMSF) with dithiothreitol (DTT)] at pH 7.4 to the ventral prostate tissues at 10 mL/g of tissue. The tissues were minced, homogenized on ice, and centrifuged for 30 min at  $30,000 \times g$  at 4°C. The supernatant was collected, pooled from all tissues, aliquoted (amounts not reported) and stored at -80°C until used. Protein concentration of the cytosol prepared for this study was determined to be 6.566 mg/mL using the Pierce BCA method (Thermo Scientific Pierce Research Lab, Rockford, IL).

**2.** <u>Saturation Radioligand Binding Experiment</u>: The summary of conditions for the saturation binding experiment is provided in Table 1 below.

TABLE 1. Summary of Conditions for Saturation Binding Experiment <sup>a</sup>				
Source of receptor		Rat prostate cytosol		
Concentration of radioligand (	(as serial dilutions)	0.25-10 nM		
Concentration of non-labeled l	igand (100X [radioligand])	25-1000 nM		
Optimization of receptor concentration		Sufficient to bind 8.6-9.0% <sup>b</sup> of radioligand at 0.25 nM		
Temperature		~2-8 °C		
Incubation time		~16 hours		
Composition of assay buffer	Tris	10 mM (pH 7.4)		
(TEDG)	EDTA	1.5 mM		
	Glycerol	10%		
	Phenylmethylsulfonyl fluoride	1.0 mM		
	DTT	1.0 mM		

a Data were not included in the study report, but are reported as a separate validation report.

On the day of the assay, the specific activity of the stock solution [ $^3$ H]-R1881 was not adjusted for decay over time, and serial dilutions in TEDG buffer were prepared to achieve the final concentrations in cytosol of 0.25, 0.5, 0.7, 1.0, 1.5, 2.5, 5.0, and 10.0 nM to determine total binding. To determine non-specific binding, solutions of non-labeled R1881 were prepared in a similar manner to achieve concentrations that were 100-fold greater than each respective radiolabeled concentration, resulting in final concentrations in cytosol of 25, 50, 70, 100, 150, 250, 500, and 1000 nM. In the absence of cytosol, the radiation found in 7.5, 15, 21, 30, or 45  $\mu$ L of 10 nM [ $^3$ H]-R1881 and 7.5, 15, or 30  $\mu$ L of 100 nM [ $^3$ H]-R1881

b As indicated in the guideline for acceptable assay performance the receptor concentration bound less than 25 to 35% of the radiolabeled R1881.

was measured. For each batch of cytosol, the optimal protein concentration was determined by calculating specific binding to differing amounts of protein per tube, using 0.25 nM radiolabeled R1881. The optimal protein concentration was determined to be 1.97 mg protein/assay tube, which resulted in the binding of 8.6-9.0% of the total radioactivity added. As indicated in the guideline for acceptable assay performance the receptor concentration bound less than 25 to 35% of the radiolabeled R1881. Cytosolic protein used in this assay was thawed fresh for this experiment at ~4°C and maintained at ~4°C during the binding assay. Each run contained three concurrent replicates at each concentration, resulting in the 72 samples depicted in Table 2.

TABLE 2. Saturation Binding Experiment Run a, b						
Total Binding	Non-Specific Binding		Radioligand alone			
Tubes 1-24 °	Tubes 25-48 d		Tubes 49-72 °			
[ <sup>3</sup> H]-R1881	[ <sup>3</sup> H]-R1881	R1881	[ <sup>3</sup> H]-R1881	[ <sup>3</sup> H]-R1881		
Final conc. (nM)	Final conc. (nM)	Final conc. (nM)	Initial conc. (nM)	$(\mu L)$		
0.25	0.25	25	10	7.5		
0.50	0.50	50	10	15		
0.70	0.70	70	10	21		
1.00	1.00	100	10	30		
1.50	1.50	150	10	45		
2.50	2.50	250	100	7.5		
5.00	5.00	500	100	15		
10.00	10.00	1000	100	30		

- a Data were not included in the study report, but are reported as a separate validation report.
- b Each concentration was run in triplicate for a total of 72 samples.
- Tubes 1-24 contained 50 μL of triamcinolone acetonide and 7.5-45 μL [3H]-R1881. Samples were dried, and 300 μL of prostate cytosol were added.
- d Tubes 25-48 contained 50 μL of triamcinolone acetonide and 7.5-45 μL [3H]-R1881. R1881 was added in a 100-fold molar excess of [3H]-R1881 in a volume of 7.5-45 μL. Samples were dried, and 300 μL of prostate cytosol were added.
- Tubes 49-72 contained only 7.5, 15, 21, 30, or 45 μL of 10 nM [3H]-R1881 or 7.5, 15, 21, or 30 μL of 100 nM [3H]-R1881 without cytosol or other components to determine the total counts added.

Following addition of triamcinolone acetonide, [3H]-R1881, and/or R1881, the tubes were dried, dissolved in diluted prostate cytosol (300 µL), and incubated for approximately 16 hours at 2-8°C. Samples were maintained at temperatures of ~4°C except during whole rack vortexing. To separate bound from free R1881, hydroxyapatite (HAP) slurry was added to each tube and vortexed once every 5 minutes for 20 minutes. The samples were then centrifuged, and the supernatant was aspirated and discarded. The samples were washed 3 times in 50 mM TRIS buffer. Following the last wash and decanting of the Tris buffer, pellets were then extracted by addition of 2 ml ethanol. The samples were vortexed 3 times at 5 minute intervals. Samples were maintained on ice at all times between vortexing. Each ethanol supernatant was then decanted into a scintillation vial, and the radiation was quantified by liquid scintillation counting. A total of 4 runs were performed on 2 batches of cytosol with similar results. For the batch of cytosol used for the competitive assay, 2 runs were performed, which had highly similar binding profiles. Final determination of acceptable AR binding assay performance was primarily based on guideline suggested standards for the competitive binding assay, although the saturation binding parameters were evaluated.

**3.** <u>Competitive Binding Experiment</u>: A summary of the assay conditions for the competitive binding experiment is included in Table 3.

TABLE 3. Summary of Conditions for Competitive Binding Experiment <sup>a</sup>					
Source of receptor		Rat ventral prostate cytosol			
Concentration of radioligand		1 nM			
Optimization of receptor conce	entration	Sufficient to bind 4.3-5.2% of 1.0 nM radioligand <sup>b</sup>			
Concentration of test substance	e (as serial dilutions)	10 <sup>-11</sup> to 10 <sup>-4</sup> M			
Incubation Temperature		4-8 °C			
Incubation time		Overnight (~16 hours)			
Composition of assay buffer	Tris	0.01 M (pH 7.4)			
	EDTA	1.5 mM			
	Glycerol	10% (v/v)			
	Phenylmethylsulfonyl fluoride with DTT	1 mM			
	Sodium molybdate	1 mM			
	Protease inhibitor	60 μΜ			

a Data were obtained from pages 15, 18, and 19 of the study report.

The competitive binding experiment was performed according to the protocol provided in the EPA Test Guidelines OCSPP 890.1150. The competitive binding experiment measures the binding of a single concentration of [<sup>3</sup>H]-R1881 (specific activity of 85.1 Ci/mmol) to the AR in the presence of increasing concentrations of a test substance. The amount of cytosolic protein used in the assay contained enough receptor to bind 4.3-5.2% of the [<sup>3</sup>H]-R1881.

Ethanol was used as the solvent vehicle, and the solubility of the test material in the vehicle and assay buffer was evaluated visually. No precipitation was noted.

Dilutions of the test substance, reference standard (R1881), weak positive control (dexamethasone), and solvent control (ethanol) were prepared to achieve the concentrations shown in Table 4. Each assay consisted of three independent runs on three different days. For each run, a set of duplicate blanks and triplicate tubes with 1 µM R1881 (non-specific binding, NSB) were run at the beginning and end of each run. Each run also included triplicate samples of each centration of the reference standard, the weak positive control, and 2,4-D, resulting in a total of 77 samples per run. In addition, duplicate blanks followed by six replicates [³H]-R1881 only (for total binding calculations) were run the day before each analysis run (the day of preparation of sample tubes).

b Data were obtained from pages 37, 39 and 41 of the study report; protein concentrations (μg/tube) were not reported.

TABLE 4. Competitor Final Molar (M) Concentrations in Competitive Binding Assay a, b					
Solvent Control	Reference standard	Weak positive control	Test Chemical	None	
Ethanol	R1881	Dexamethasone	2,4-D	None	
Tubes 3-5 and 72-74	Tubes 6-23 and 75-77 °	Tubes 24-47	Tubes 48-71	Tubes 1-2	
	1×10 <sup>-6</sup>	$1 \times 10^{-3}$	$1 \times 10^{-4}$		
	$1 \times 10^{-7}$	$1 \times 10^{-4}$	$1 \times 10^{-5}$		
	$1 \times 10^{-8}$	1×10 <sup>-5</sup>	$1 \times 10^{-6}$		
	1×10 <sup>-9</sup>	1×10 <sup>-6</sup>	$1 \times 10^{-7}$		
	$1 \times 10^{-10}$	$1 \times 10^{-7}$	$1 \times 10^{-8}$		
	1×10 <sup>-11</sup>	1×10 <sup>-8</sup>	1×10 <sup>-9</sup>		
		1×10 <sup>-9</sup>	$1 \times 10^{-10}$		
		$1 \times 10^{-10}$	$1 \times 10^{-11}$		

- a Data were obtained from pages 37-42 of the study report.
- b Each concentration of each chemical was run in triplicate, plus duplicate blanks for a total of 77 tubes per run. Tubes 3-77 contained 50 μL of triamcinolone acetonide and 30 μL [³H]-R1881. Samples were dried, and 300 μL of prostate cytosol were added. Tubes 3-77 also contained 10 μL of the solvent control, reference standard (non-radiolabeled R1881), weak positive control, or test substance, with the exception of Tubes 6-8 and 75-77 that contained 30 μL of non-radiolabeled R1881 (used to evaluate non-specific binding). Six tubes analyzed the day prior to each run analysis contained only 30 μL of [³H]-R1881 to determine ligand activity.
- c Tubes 6-8 and 75-77 were used to evaluate non-specific binding by adding 100x of cold (non-radiolabeled) R1881.

Sample tubes were stored overnight at 4-8°C in the dark to allow the reaction to reach equilibrium, bound R1881 was separated from free R1881 by washing with HAP buffer and extraction with ethanol, followed by scintillation counting of bound [<sup>3</sup>H]-R1881.

**4.** <u>Data Analysis</u>: The top and bottom of the curve, Hill slope, inhibition concentration (IC<sub>50</sub>), and standard deviations were assessed using GraphPad Prism v. 5, and the data were fitted to a "one site binding" non-linear regression model (GraphPad Prism v. 5).

# 5. Definitions

# a. Classification of test material

If the data fit a 4-parameter nonlinear regression model, the test chemical is classified as:

**Binder:** The average curve for the test chemical across runs crosses 50% of radioligand bound.

**Equivocal:** The average lowest portion of curves across runs is between 50% and 75% radioligand binding (*i.e.* radioligand displacement is at least 25% but less than 50%), or the curve falls outside the range for the weak positive control (-0.6 to -1.4).

**Non-Binder:** The average lowest portion of curves across runs is greater than 75% activity (*i.e.* less than 25% displacement of radioligand), or the data do not fit the model.

**Untestable:** If the test compound is not soluble above  $1 \times 10^{-6}$  M and the binding curve does not cross 50%, the chemical is judged to be untestable.

Table 5. Data Interpretation Criteria Specified in OCSPP 890.1150					
	Classification	Values			
	Average curve across runs crosses 50% <sup>a</sup>	Binder	2		
Data fit 4-parameter nonlinear regression	Average lowest portion of curves across runs is between 50% and 75% activity <sup>b</sup>	Equivocal	1		
model	Average lowest portion of curves across runs is greater than 75% activity <sup>b</sup>	Non-Rinder	0		
Data do not fit the model		Non-Binder			

- a If the curve fell outside the range for the weak positive control (see test acceptability criteria), the run as classified as equivocal.
- b If the test compound was not soluble above 10<sup>-6</sup> M and the binding curve did not cross 50%, the chemical was judged to be untestable.

# b. Descriptors for receptor binding

Bmax: maximal binding capacity

**K**<sub>d</sub>: dissociation constant

IC50: Concentration of the test substance at which 50% of radioligand is displaced from the AR

by the competitor

**Relative Binding Affinity (RBA):**  $IC_{50}$  of R1881 × 100 ÷  $IC_{50}$  of test substance

#### II. RESULTS

**A.** <u>SATURATION BINDING EXPERIMENT</u>: Saturation binding experiment parameters are presented in Table 6. The dissociation constant ( $K_d$ ) for [ $^3H$ ]-R1881 was 0.4641, and the estimated  $B_{max}$  (nM) was 0.06392 for the batch of prostate cytosol that was used for this study. The  $K_d$  was below the recommended range reported in the test guideline (0.685-1.57 nM). Confidence in these numbers is high according to the goodness of fit ( $R^2 = 0.9871-0.9937$ ) and the small variation among runs.

TABLE 6. Saturation Binding Experiment of R-1881 with Androgen Receptor from Rat Prostate						
Cytosola						
Parameter	Run 1 <sup>b</sup>	Run 2 <sup>b</sup>	Run 3 b	Mean Runs 1-2 °		
R <sup>2</sup> (unweighted)	0.9937	0.9871	ND	0.9871-0.9937		
$B_{max}$ (nM)	0.06194	0.06590	ND	0.06392		
B <sub>max</sub> (fmol/100μg protein)	3.146	3.343	ND	3.245		
$K_{d}(nM)$	0.4359	0.4922	ND	0.4641		

- Data were not included in the study report, but are reported as a separate validation report.
- b Two saturation runs were performed for this batch of cytosol.
- c The range of  $R^2$  is reported and the mean is reported for the other parameters.
- R<sup>2</sup> Goodness of fit for curve calculated for specific binding,
- ND Not determined

Figure 1 illustrates the non-specific, specific, and total binding curves for [<sup>3</sup>H]-R1881 to the androgen receptor. The specific binding reached a plateau and the non-specific binding was generally less than 20% of total binding at all concentrations (range 6.2%-19.8%) except the highest concentration in Run 1 (24.6%). All other values indicated acceptable performance of

the assay. Figure 2 is a Scatchard plot that illustrates the binding of [³H]-R1881 to the androgen receptor. The data fit results in a linear plot.

FIGURE 1. Binding of [<sup>3</sup>H]-R1881 to the Androgen Receptor during the Saturation Binding Experiment.


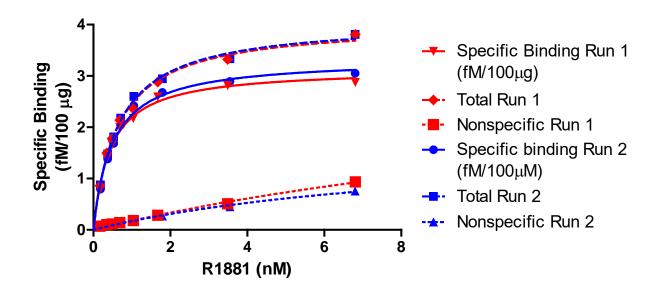


FIGURE 2. Scatchard Plot of the Binding of [3H]-R1881 to the Androgen Receptor.

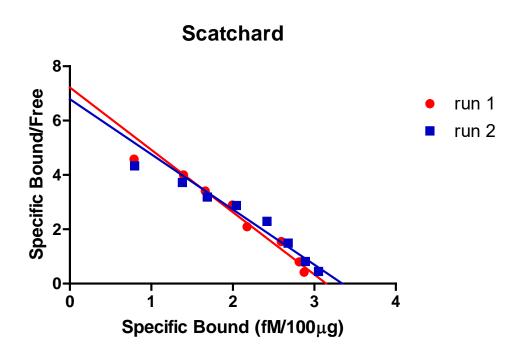

**B.** COMPETITIVE BINDING EXPERIMENT: The results from the three competitive binding experiments are summarized in Tables 7 and 8 and shown graphically in Figures 3 and 4. The estimated mean log IC<sub>50</sub> for 2,4-D could not be determined as it did not result in 50% displacement of the radioligand at any concentration. The estimated average log IC<sub>50</sub>s for R1881 and the weak positive control (dexamethasone) were −9.0 and −4.4 M, respectively. The mean RBA for the positive control was 0.0027%. Confidence in these numbers is high due to the small variation. No precipitation of the test compound was visually observed at any concentration (≤10<sup>-4</sup> M). The solvent control responses indicated no drift in the study assay.

TABLE 7. Competitive Binding Experiment Results for Strong and Weak Positive with AR from Rat Prostate Cytosol <sup>a,b</sup>

Prostate Cytosor	Assay #1	Assay #2	Assay #3	Mean ± SEM	Performance Criteria
R1881		it (%) of Total Bi	nding		
10-11	99.8	96.3	101.7	99.0±1.6	
10 <sup>-10</sup>	92.1	87.3	93.5	90.9 ±1.9	
10-9	50.0	51.1	52.5	51.2 ±0.7	
10-8	10.0	9.5	9.7	9.7 ±0.2	
10-7	0.2	0.2	0.1	0.2 ±0.0	
Bottom (%)	-0.6	-1.6	-1.1	-1.1	-2.0 to 2.0
Top (%)	101	97	103	101	82 to 114
Log IC <sub>50</sub> (M)	-8.993	-8.986	-8.956	-8.978	
IC <sub>50</sub> (M)	1.017e-009	1.032e-009	1.106e-009	1.052e-009	
Hill slope	-1.0	-0.9	-1.0	-1.0	-1.2 to -0.8
$R^2$	0.9998	0.9999	1.000	1.000	
RBA (%)	100	100	100	100	100
Dexamethasone	Percen	t (%) of Total Bi	nding		
10 <sup>-10</sup>	98.2	97.9	104.2	100.1±2.0	
10 <sup>-9</sup>	98.1	97.5	104.0	99.9 ±2.1	
10-8	97.9	99.0	102.1	99.7 ±1.3	
10-7	96.5	98.0	104.7	99.7±2.5	
10 <sup>-6</sup>	94.7	96.3	99.3	$96.8 \pm 1.3$	
10 <sup>-5</sup>	79.6	79.1	80.5	$79.7 \pm 0.4$	
10 <sup>-4</sup>	29.7	25.3	27.7	$27.6 \pm 1.3$	
10 <sup>-3</sup>	3.9	2.3	3.6	3.2 ±0.5	
Bottom (%)	-0.5	-0.5	-0.6	-0.6	-12 to 12
Top (%)	98	98	104	100	87 to 106
Log IC <sub>50</sub> (M)	-4.379	-4.440	-4.416	-4.412	
IC <sub>50</sub> (M)	4.177e-005	3.630e-005	3.840e-005	3.871e-005	
Hill slope	-1.0	-1.1	-1.0	-1.0	-1.4 to -0.6
$R^2$	0.9998	0.9999	0.9995	1.000	
RBA (%)	0.0024	0.0028	0.0029	0.0027	

a Data were obtained from page 30 of the study report.

b The mean and standard deviations are reported for the combined runs.

NA Not applicable. R<sup>2</sup> Goodness of fit

IC<sub>50</sub> Concentration of the test substance at which 50% of radioligand is displaced from the AR by the competitor

RBA (%) Relative binding affinity

TABLE 8. Competitive Binding Experiment Results for 2,4-D with AR from Rat Prostate Cytosol a,b
---

TABLE 6. Competiti	Assay #1	Assay #2	Assay #3	Mean ± SEM
2,4-D	Perc	ent (%) of Total Bi	nding	
10-11	98.79	98.24	99.89	$99.0 \pm 0.5$
$10^{-10}$	102.44	99.72	99.55	$100.6 \pm 0.9$
10-9	103.53	97.07	99.95	$100.2 \pm 1.9$
10-8	103.22	99.22	104.25	$102.2 \pm 1.5$
10 <sup>-7</sup>	101.49	99.04	103.98	$101.5 \pm 1.4$
10 <sup>-6</sup>	103.17	98.32	104.51	$102.0 \pm 1.9$
10 <sup>-5</sup>	102.33	97.65	105.58	$101.9 \pm 2.3$
$10^{-4}$	105.64	99.8	106.32	$103.9 \pm 2.1$
Bottom (%)	103.2	n/a	99.70	n/a
(95% CI)	(101.1 to 105.4)		(97.9 to 101.5)	
Top (%)	~-146.6	n/a	105.1	n/a
(95% CI)	(very wide)		(103.8 to 106.4)	
$Log IC_{50}(M)$	n/a	n/a	n/a	n/a
(95% CI)				
$IC_{50}(M)$	n/a	n/a	n/a	n/a
(95% CI)				
$Log EC_{50}(M)$	~-13.15	n/a	-8.363	n/a
(95% CI)	(very wide)		(-9.431 to -7.294)	
$EC_{50}(M)$	~7.139e-014	n/a	4.339e-009	n/a
(95% CI)	(very wide)		(3.707e-010 to	
			5.078e-008)	
Hill slope	~-0.8113	n/a	1.939	n/a
(95% CI)	(very wide)		(-3.018 to 6.895)	
$R^2$	0.6318	n/a	0.9371	n/a
Data Interpretation	negative	negative	negative	n/a
RBA (%)	n/a	n/a	n/a	n/a

a Data were obtained from page 31 of the study report.

b The mean and standard deviations are reported for the combined runs.

n/a Not applicable. R<sup>2</sup> Goodness of fit

IC<sub>50</sub> Concentration of the test substance at which 50% of radioligand is displaced from the AR by the competitor

RBA (%) Relative binding affinity

## **Collective Responses of the Independent AR Binding Assays**

The collective responses of the three independent AR binding assays for 2,4-D indicates no apparent alterations in radiolabeled R1881 binding at any of the concentrations tested. The mean results of the three independent assays are shown in Figure 4 (R1881 and dexamethasone controls) and Figure 5 (2,4-D) of the study report. The final classification of 2,4-D was based on the average of the three valid assays with a mean value of "0" (i.e., (0+0+0)/3) as described in Text Table 5 of the study report. Based on these data, 2,4-D was classified as non-binding at concentrations up to  $10^{-4}$  M.


## **Data Interpretation**


## Test Validity

The QC criteria for the reference chemicals in the AR binding assay indicated that the assay performed according to the specified criteria. As with any biological system, there was slight

variability between assays, but the overall robustness of the responses for R1881 (strong positive control/standard curve) and dexamethasone (weak positive control) indicated that each assay included in this assessment performed as expected. Thus, assay #1, #2, and #3 of the AR binding assay with 2,4-D were considered valid.

FIGURE 3. Percentage R1881 Bound to the Androgen Receptor in the Presence of Radioinert R1881, Dexamethasone, and 2,4-D (Assays 1 – 3).

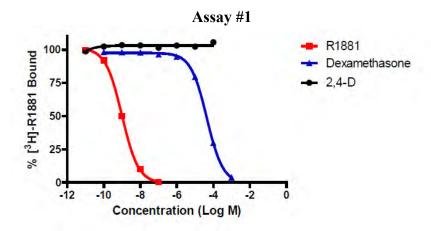

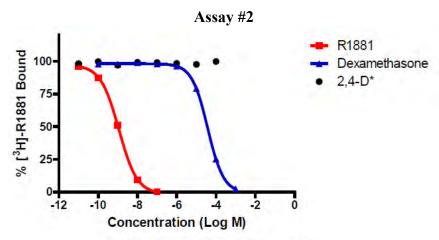

\*Non-linear regression could not be calculated.


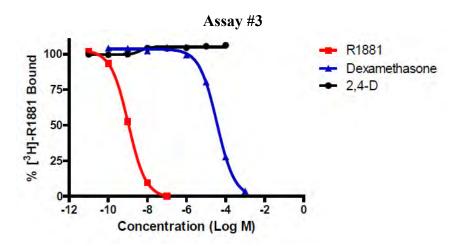


FIGURE 4. Mean of Percentage R1881 Bound to the Androgen Receptor in the Presence of 2,4-D from Three Assays.

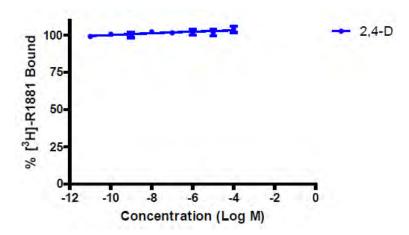

C. PERFORMANCE CRITERIA: To ensure that the competitive binding assay was functioning properly, each run was evaluated using the criteria in Table 9. The curve for the reference material showed that increasing concentrations of unlabeled R1881 displaced [<sup>3</sup>H]-R1881 in a manner consistent with one-site binding, as indicated by a Hill slope of -1.0 to -0.9. Examination across the runs indicated consistency of the Hill slope, placement along the X-axis, and top and bottom plateaus.

TABLE 9. Criterion <sup>a</sup>	Tolerance Limit(s) b	Value	Yes	No
<b>Ligand depletion</b> is minimal. The recommended ratio of				
total binding in the absence of competitor to total amount	≤15%	4.3-5.2%	X	
of [3H]-R1881 added per assay tube.				
<b>2,4-D</b> Top (% binding)	80 to 115	100 to 106	X	
R1881 fitted curve parameters				
Top (% binding)	82 to 114	97 to 103	X	
Bottom (% binding)	-2.0 to 2.0	−1.6 to −0.6	X	
Hill Slope	−1.2 to −0.8	−1.0 to −0.9	X	
Weak positive control (dexamethasone) fitted curve param	neters			
Top (% binding)	87 to 106	98 to 104	X	
Bottom (% binding)	-12 to 12	−0.6 to −0.5	X	
Hill Slope	−1.4 to −0.6	-1.1 to -1.0	X	
Saturation Binding Experiment Kd (nM)	(0.685-1.57 nM)	0.4641		X
Non-specific binding (%)	≤10.0	7.30 to 8.23 °	X	

a Data were obtained from pages 30, 31 and 37-42 of the study report.

NR Not reported

b These values represent ranges from the validation study.

c Values reported for the three NSB tubes at the beginning of each run; does not include the three NSB tubes at the end of the run.

#### III. DISCUSSION AND CONCLUSIONS

- A. <u>INVESTIGATOR'S CONCLUSIONS</u>: Based on the combined responses in each of three independent androgen receptor binding assays, it was determined that 2,4-D had no appreciable effect in the binding of the reference androgen ([3H]-R1881) at any concentration (up to 10<sup>-4</sup> M). The results of the *in vitro* AR binding assay using rat prostate cytosol indicate that, under the conditions of this study, 2,4-D was negative for AR binding at concentrations up to 10<sup>-4</sup> M.
- **B.** AGENCY COMMENTS: The saturation binding experiment resulted in a maximum binding capacity (B<sub>max</sub>) of 3.245 fmol/100 μg protein and the dissociation constant (K<sub>d</sub>) was 0.4641 nM. Although these values were slightly below the range of values from the validation studies, the results were highly reproducible and all other performance criteria and the competitive binding assays indicated acceptable performance of the assay

The test guideline recommends testing up to  $10^{-3}$  M; however, the sponsor selected  $10^{-4}$  M as the highest concentration, based on *in vivo* toxicokinetic analyses in the rat. Specific binding was >75% at all concentrations tested ( $10^{-11}$  to  $10^{-4}$  M). An IC<sub>50</sub> and RBA could not be calculated for 2,4-D as it did not result in 50% displacement at any concentration.

The estimated average log IC<sub>50</sub>s for R1881 and the weak positive control (dexamethasone) were -9.0 and -4.4 M, respectively. The mean RBA for the positive control was 0.0027%. Confidence in these numbers is high due to the small variation. No precipitation of the test compound was visually observed at any concentration ( $\le 10^{-4}$  M). The solvent control responses indicated no drift in the study assay, and all performance criteria were met in all three runs.

Based on the results of the three runs, 2,4-D is classified as a non-binder for the androgen receptor at concentrations up to  $10^{-4}$  M.

- C. <u>STUDY DEFICIENCIES</u>: The following deficiencies were noted that are not considered to have had an adverse impact on the results, interpretation or conclusions of this study:
  - Only two saturation binding runs were conducted rather than the three runs recommended in the test guideline.

# DATA EVALUATION RECORD

2,4-DICHLOROPHENOXYACETIC ACID (2,4-D)

Study Type: OCSPP 890.1200, Aromatase Assay

EPA Contract No. EP10H001452 Task Assignment No. 2-26-2012 (MRID 48614302)

Prepared for
Health Effects Division
Office of Pesticide Programs
U.S. Environmental Protection Agency
2777 South Crystal Drive
Arlington, VA 22202

Prepared by CSS-Dynamac Corporation 1910 Sedwick Road, Building 100, Suite B Durham, NC 27713

Primary Reviewer:	Signature:	David a. M. Ewer
David A. McEwen, B.S.	Date:	3/17/2012
Secondary Reviewer:	Signature:	Lucientes
Scott D. Studenberg, Ph.D., D.A.B.T.	Date:	3/23/2012
D 14	G:	Jack Q. Enry
Program Manager:	Signature:	
Jack D. Early, M.S.	Date:	3/26/2012
		Jack Q. Ewy
Quality Assurance:	Signature:	
Jack D. Early, M.S.	Date:	3/26/2012

This Data Evaluation Record may have been altered by the Health Effects Division subsequent to signing by CSS-Dynamac Corporation personnel.

The US EPA Endocrine Disruptor Screening Program (EDSP) Tier 1 screening battery is comprised of eleven screening assays intended to identify a chemical's likely endocrine bioactivity, i.e., its potential to interact with the estrogen, androgen, or thyroid (E, A, or T) pathways. The robustness of the Tier 1 battery is based on the strengths of each individual assay to identify potential endocrine bioactivity with complementary endpoints within the assay, where available, and redundancy across the battery. Thus, the results of each individual assay should not be considered in isolation but rather should be considered in the context of other assays in the battery as well as Other Scientifically Relevant Information (OSRI). In order to determine if a chemical has the potential to interact with the E, A or T pathways, a Weight of Evidence (WoE) evaluation of Tier 1 assay results, in combination with the findings in the OSRI, should be undertaken (refer to the WoE Document).

Primary Reviewer: Patience Browne, Ph.D.

Office of Science Coordination and Policy

Secondary Reviewer: Greg Akerman, Ph.D.

Health Effects Division, Office of Pesticide Programs

Signature: Z

Signature:

Date: 6/9/15
Template version 08/2011

# DATA EVALUATION RECORD

STUDY TYPE: Aromatase (Human Recombinant); OCSPP 890.1200

PC CODE: 030001 DP BARCODE: D398640

TXR#: 0052104 CAS No.: 94-75-7

TEST MATERIAL (PURITY): 2,4-D (98.5% a.i.)

**SYNONYMS**: 2,4-Dichlorophenoxyacetic acid

**CITATION:** Coady, K.K. and Sosinski, L.K. (2011) 2,4-Dichlorophenoxyacetic acid:

Evaluation of 2,4-dichlorophenoxyacetic acid in the human recombinant aromatase assay. Toxicology & Environmental Research and Consulting, The Dow Chemical Company, Midland, MI. Laboratory Project Study ID: 111036,

September 27, 2011. MRID 48614302. Unpublished.

SPONSOR: Industry Task Force II on 2,4-D Research Data, c/o McKenna Long & Aldridge

LLP, 1900 K Street NW, Washington, D.C.

TEST ORDER #: CON-030001-01

EXECUTIVE SUMMARY: In an *in vitro* aromatase (CYP 19) assay (MRID 48614302), 2,4-D (98.5% a.i., lot#: 2006 2433 8006-USA) in ethanol (1%) was incubated with human recombinant aromatase and tritiated androstenedione (1-β [³H(N)]-androst-4-ene-3,17-dione; [³H]ASDN) for 15 minutes at 37 °C to assess the potential of 2,4-D to inhibit aromatase activity. 2,4-D was tested at logarithmic concentrations from 10<sup>-10</sup> M to 10<sup>-4</sup> M in three independent runs. The test guideline recommends testing up to 10<sup>-3</sup> M; however, the sponsor selected 10<sup>-4</sup> M as the highest concentration based on *in vivo* toxicokinetic analyses in the rat. The sponsor considered concentrations higher than 10<sup>-4</sup> M to be not relevant for testing in this assay as they are substantially above the inflection point for linear toxicokinetics (See Appendix B of the AR binding assay study report, MRID 48614301).

Aromatase activity was determined by measuring the amount of tritiated water produced at the end of a 15-minute incubation for each concentration of chemical. Tritiated water was quantified using liquid scintillation counting (LSC). Three independent runs were conducted and each run included a full activity control, a background activity control, a positive control series ( $10^{-10}$  to  $10^{-5}$  M) using a known inhibitor (4-hydroxyandrostenedione; 4-OH ASDN), and the test chemical series ( $10^{-10}$  to  $10^{-4}$  M) with 3 repetitions per concentration.

Aromatase activity in the full activity controls ranged from 0.131 to 0.244 nmol·mg-protein<sup>-1</sup>·min<sup>-1</sup> for the three test runs, with a mean and standard deviation of  $0.186\pm0.036$  nmol·mg-protein<sup>-1</sup>·min<sup>-1</sup>. Activity in the background controls ranged from 7.31 to 11.51% of the full activity controls. The responses of full activity controls were outside of the 90 to 110% range in the  $2^{nd}$  and  $3^{rd}$  replicates in the Run 1 (86.6 and 113.8%, respectively).

Results for the positive control were generally within the recommended ranges for the top of the curve, bottom curve, Hill slope, log IC<sub>50</sub>, and coefficient of variation for replicates of each concentration within runs, with the exception that the bottom of the curve in Runs 2 and 3 (-8.2 and -7.2, respectively) exceeded the acceptable range (-5 to +6). Also, the coefficients of variations (%CVs) for replicates of each concentration of 2,4-D within a run were generally within the 15% guideline, with the exception that the overall %CV for the highest two concentrations (-45.2 and 415.9%) exceeded the acceptable limit. For 4-OH ASDN, the estimated log IC<sub>50</sub> averaged -7.23 M and the slope was -0.92.

For 2,4-D, aromatase activity averaged 0.195±0.045 nmol·mg-protein<sup>-1</sup>·min<sup>-1</sup> at the lowest tested concentration (10<sup>-10</sup> M) and 0.168±0.025 nmol·mg-protein<sup>-1</sup>·min<sup>-1</sup> at the highest tested concentration (10<sup>-4</sup> M). The data for 2,4-D were modeled for Run 1, and the goodness of fit (R<sup>2</sup>) value was 0.83; however, the 2,4-D data from Runs 2 and 3 could not be modeled. The average dose-response curve indicated that the aromatase activity of the test material at concentrations ranging from 10<sup>-10</sup> M to 10<sup>-4</sup> M was essentially equivalent to the activity observed in the full activity controls. At 10<sup>-4</sup> M, aromatase activity was approximately 91%. Since the average lowest portion of the activity response curve was greater than 75% activity, 2,4-D is classified as a non-inhibitor of aromatase activity up to the highest concentration tested (10<sup>-4</sup> M). High CVs were observed for 4-OH ASDN at the two highest concentrations and at a single concentration for 2,4-D. Individual values were occasionally outside of the performance criteria ranges (with the mean value within range).

Based on the data from the average response curve, 2,4-D is classified as a Non-inhibitor of aromatase activity in this assay.

The assay **satisfies** the EDSP Tier 1 Test Order requirements for an Aromatase assay (OCSPP 890.1200).

**<u>COMPLIANCE</u>**: Signed and dated Data Confidentiality, GLP Compliance, and Quality Assurance statements were provided.

#### I. MATERIALS AND METHODS

## A. MATERIALS

**Structure:** 

1. Test Substance: 2,4-D

**Description:** Off-white powder

Source: Nufarm Americas, Inc. (Burr Ridge, IL)

**Lot/Batch #:** 2006 2433 8006-USA

Purity: 98.5% a.i.
Volatility: Not reported
Storage conditions: Ambient
Stability: Not reported
Solvent: Ethanol

**Solubility (in test solvent):** Soluble up to 0.01 M **Highest Concentration Tested:**  $10^{-4} \text{ M in all runs}$ 

**Stock Solution Preparation** 

**Methodology:** Dissolved the test material in ethanol.

**Molecular weight:** 221.0 **CAS #:** 94-75-7

OH CI CI

2. Non-Labeled Substrate: Androstenedione (ASDN)

CAS #: 63-05-8

Source: Steraloids, Inc. (Cat. # A6030100)

**Lot/Batch #:** L1627 **Purity:** 98.4%

3. Radiolabeled Substrate:  $1-\beta$  [ $^3$ H(N)]-Androst-4-ene-3,17-dione; ([ $^3$ H]ASDN)

Source: Perkin Elmer Life and Analytical Sciences (Cat. # NET 926)

Lot/Batch #: 619344
Radiochemical Purity (Supplier): >97%
Specific activity: 26.3 Ci/mmol

Radiochemical Purity (In-lab

determination): Not reported

**4. Positive Control:** 4-hydroxyandrostenedione (4-OH ASDN)

**CAS #:** 566-48-3

Source: Sigma-Aldrich (Cat. # F2552)

**Lot/Batch #:** 081K2133 **Purity:** 99.6%

**5. Solvent (Vehicle Control):** Ethanol

Source: Sigma Aldrich (Cat. # E7023)

Lot/Batch #: 72596PMV

**Justification for choice of solvent:** Not provided. Ethanol is a Guideline preferred solvent.

Concentration

(% of total volume in assays): 1% v/v

6. Test Microsomes: Human recombinant aromatase (CYP19) microsomes

Source: Gentest (Woburn, MA; Cat. # 456260)

Lot/Batch #:03897Protein concentration:7.4 mg/mLCytochrome C reductase activity:290 nmol/min/mgAromatase activity:6.0 pmol/min/pmol P450

# B. METHODS

1. <u>Assay Components and Preparations</u>: A mixture of non-labeled and radiolabeled [<sup>3</sup>H]-ASDN was prepared to result in a 2 μM ASDN solution with a predicted radioactive content of 1.0 μCi/mL.

Test chemical stock solutions were prepared such that the total volume of each test chemical formulation used per assay was no more than 1% v/v of the total assay volume. The report specified that ethanol was chosen because it was mentioned in the guideline as a preferred solvent.

A stock solution of the positive control substance, 4-OH ASDN, was formulated in ethanol. Fresh serial dilutions of the stock solution were prepared each time the aromatase inhibition assay was conducted. Dilutions were prepared such that the target concentrations of the positive control substance ( $10^{-10}$  to  $10^{-5}$  M; Table 4) were achieved by the addition of 20  $\mu$ L of the dilution for a final assay volume of 2 mL.

Human recombinant microsomes were purchased from Gentest, and aliquoted into individual vials based on protein content. Microsomes were stored at approximately -80° C until use.

Other assay components sodium phosphate buffer, propylene glycol, and NADPH are reported in Table 1.

TABLE 1. Assay Components and Conditions				
Assay Factor	Values			
0.1M sodium phosphate buffer (pH 7.4)				
Microsomal Protein	0.004 mg/mL			
NADPH	6 mM			
[³H]ASDN	2 μΜ			
Propylene Glycol	100 μL			
Temperature	37°C			
Incubation Time	15 min			

- 2. <u>Suitability Assessments</u>: The protein concentration in an aliquot of the microsomes was determined each day of use, and microsomes were diluted with phosphate buffer such that approximately 0.004 mg/mL protein was present in the final reaction solution. Aromatase activity of the microsomes was provided by the vendor as 6.0 pmol/min/pmol P450. The minimum aromatase activity in the full activity control samples was determined to be 0.131 nmol/mg-protein/min, which was greater than the minimum acceptable aromatase activity of 0.10 nmol/mg-protein/min.
- **3.** Aromatase Assay: Each assay run contained four tubes for the full enzyme activity and four tubes for the background activity controls. Two tubes of each control were run at the

beginning of the assay, and two of each control were run at the end of the assay. A full concentration curve in duplicate for the positive control, and a full concentration curve in triplicate for the test substance were established. The aromatase assay was conducted according to the procedures described in OCSPP 890.1200 (Section h, pp. 9-10).

The amount of  ${}^3H_2O$  in the aqueous fraction was quantified for each assay tube by LSC, and aromatase activity was reported in units of nmol·mg-protein ${}^{-1}$ ·min ${}^{-1}$ .

**Demonstration of Proficiency:** It was stated that all assays were performed by personnel with demonstrated proficiency performing the assay as outlined in test guideline 890.1200. Proficiency records are archived in training records at the test facility and a copy of these data were archived with the 2,4-D specific report.

## a. **Positive Control**

- (1) <u>Initial Demonstration of Laboratory Proficiency</u>: Raw data were not provided, however, summary data of proficiency exercised were equivocal. Relative aromatase activity in the presence of prochloraz, a known aromatase inhibitor indicated dose-dependent inhibition. Aromatase activity appeared to be unaffected by ronidazole, a non-inhibitor at concentrations ranging from 10<sup>-10</sup> to 10<sup>-2.5</sup> M. Dose-dependent inhibition of aromatase was demonstrated in the presence of fenarimol and nitrofen, two recognized inhibitors of aromatase; however, the data were highly variable.
- **Demonstration of Proficiency of New Technician for Conducting Assay (when applicable):** Demonstration of proficiency by a new technician, if applicable, was not reported. The positive control data for slope, top and bottom percent, and log IC<sub>50</sub> met the criteria as listed in section (i) of OCSPP 890.1200, with the exception of the bottom of the curve, which was below the recommended value for Run 2 (–8.2) and Run 3 (–7.2).

TABLE 2. Performance Criteria for the Positive Control						
Parameter	Lower Limit Criteria	Upper Limit Criteria	Actual Lower Limita	Actual Upper Limita		
Slope	-1.2	-0.8	-1.1	-0.78		
Top (%)	90	110	91	100		
Bottom (%)	-5	+6	-8.2	-2.6		
Log IC <sub>50</sub> (M)	-7.3	-7.0	-7.3	-7.2		

a Data were obtained from page 32 of the study report.

**b.** <u>Proficiency Chemicals</u>: Standard curves were provided for the reference chemicals prochloraz, ronidazole, fenarimol, and nitrofen (Appendix A of the protocol on pages 84-86 of the study report).

TABLE 3. Proficiency Chemicals					
Compound	CAS#	Class	Concentrations <sup>a</sup>		
Prochloraz	Not provided	Inhibitor	10 <sup>-10</sup> to 10 <sup>-3</sup> M		
Fenarimol	60168-88-9	Inhibitor	10 <sup>-10</sup> to 10 <sup>-3</sup> M		
Nitrofen	1836-75-5	Inhibitor	10 <sup>-10</sup> to 10 <sup>-3</sup> M		
Ronidazole	Not provided	Non-inhibitor	10 <sup>-10</sup> to 10 <sup>-3</sup> M		

a Concentration ranges taken from figures.

5. Determination of Aromatase Activity with Test Chemical(s): The response of aromatase activity to the presence of eight concentrations of 2,4-D per run, in triplicate, was tested during three independent runs (Table 4). Solubility was visually assessed (presence of cloudiness or a precipitate). No precipitation was observed at any concentration in any run. The response of each full activity control within a run was between 91 to 108% of the average full activity, with the exception of two full activity control responses (87 and 114%) from the beginning and end, respectively, of Run 1 that were outside of the guideline recommended range of 90 to 110%.

TABLE 4. Test Chemical St	ABLE 4. Test Chemical Study Design for each Test Run				
Sample Type	Repetitions (Tubes)	Description	Reference or Chemical (M)		
Full Activity Control	4	All test components <sup>a</sup> plus solvent vehicle	N/A		
Bkgd Activity Control	4	Same as above without NADPH	N/A		
4-OH ASDN Conc 1	2	All test components plus 4-OH ASDN	1×10 <sup>-5</sup>		
4-OH ASDN Conc 2	2	All test components plus 4-OH ASDN	1×10 <sup>-6</sup>		
4-OH ASDN Conc 3	2	All test components plus 4-OH ASDN	1×10 <sup>-6.5</sup>		
4-OH ASDN Conc 4	2	All test components plus 4-OH ASDN	1×10 <sup>-7</sup>		
4-OH ASDN Conc 5	2	All test components plus 4-OH ASDN	1×10 <sup>-7.5</sup>		
4-OH ASDN Conc 6	2	All test components plus 4-OH ASDN	1×10 <sup>-8</sup>		
4-OH ASDN Conc 7	2	All test components plus 4-OH ASDN	1×10 <sup>-9</sup>		
4-OH ASDN Conc 8	2	All test components plus 4-OH ASDN	1×10 <sup>-10</sup>		
2,4-D Conc 1	3	All test components plus 2,4-D	1×10 <sup>-4</sup>		
2,4-D Conc 2	3	All test components plus 2,4-D	1×10 <sup>-4.5</sup>		
2,4-D Conc 3	3	All test components plus 2,4-D	1×10 <sup>-5</sup>		
2,4-D Conc 4	3	All test components plus 2,4-D	1×10 <sup>-6</sup>		
2,4-D Conc 5	3	All test components plus 2,4-D	1×10 <sup>-7</sup>		
2,4-D Conc 6	3	All test components plus 2,4-D	1×10 <sup>-8</sup>		
2,4-D Conc 7	3	All test components plus 2,4-D	1×10 <sup>-9</sup>		
2,4-D Conc 8	3	All test components plus 2,4-D	1×10 <sup>-10</sup>		

a The complete assay contained buffer, propylene glycol, microsomal protein, Γ<sup>3</sup>H]ASDN, and NADPH.

# C. <u>DATA ANALYSIS</u>

1. <u>Raw Data:</u> Raw data were converted to aromatase activity (nmol/mg protein/min) and percent control for the positive control and test chemical. The following raw data and calculated endpoints for each run were included in the report (Table 5).

TABLE 5. Raw and Calculated Data	
Raw/Calculated Data	Included (X)
DPM/mL for each portion of extracted aqueous incubation mixture	X
Average DPM/mL for each aqueous portion (after extraction)	X
Total DPM for each aqueous portion (after extraction)	X
The total DPM present in the assay tube at initiation	X
The percentage of substrate converted to product	X
Total DPM after extraction corrected for background	X
Aromatase activity expressed in nmol/mg protein/min	X
Average aromatase activity in the full activity control tubes	X
Percentage of control activity remaining in the presence of various inhibitor concentrations	X

DPM Disintegrations per minute

**Statistical Methods:** Statistical analyses and graphical displays were conducted using Graph Pad Prism (Version 4.0, La Jolla, CA). Basic statistical analyses were performed on the data, which included means of replicates, standard deviation of the mean, relative standard deviation, and coefficient of variation. The Hill slope and log IC<sub>50</sub> values across three independent runs were compared based on a one-way random effects analysis of variance, treating runs as random effects.

The response curve was fitted by nonlinear regression analysis. Model fits were carried out using a 4-parameter regression model. For each run, percent of full activity control were plotted versus logarithm (base 10) of the test chemical concentration or 4-OH ASDN concentration. Each run was plotted with the data's best fit curve. Additionally, the average inhibition response curve across all runs was also plotted.

**3.** <u>Interpretation of Results</u>: Interpretation of the assay results was based on the average of three runs, using the categories presented in Table 6.

TABLE 6. Interpretation of Results					
	Criteria	Interpretation			
Data fit 4-parameter nonlinear	Average curve across runs crossed 50% <sup>a</sup>	Inhibitor			
regression model	Average lowest portion of curves across runs is between 50% and 75% activity <sup>b</sup>	Equivocal			
	Average lowest portion of curves across runs is greater than 75% activity <sup>b</sup>	Non-inhibitor			
Data do not fit model					

a Ordinarily, an inhibition curve will fall from 90% to 10% over 2 log units with a slope near -1. Unusually steep curves may indicate protein denaturing or solubility issues. If the slope of the curve is steeper than -2.0, the result is classified as equivocal.

#### II. RESULTS

- A. <u>CONTROL ACTIVITY</u>: Aromatase activity in the full activity controls ranged from 0.131 to 0.244 nmol·mg-protein<sup>-1</sup>·min<sup>-1</sup> for the three test runs, with a mean and standard deviation of 0.186±0.036 nmol·mg-protein<sup>-1</sup>·min<sup>-1</sup>. Activity in the background controls ranged from 7.31 to 11.51% of the full activity controls. The response of each full activity control was generally between 90 to 110% of the average full activity, with the exception of the 2<sup>nd</sup> and 3<sup>rd</sup> replicates in the Run 1 (86.6 and 113.8%, respectively). The response of the full activity controls and background controls were acceptable.
- **B. POSITIVE CONTROL:** For the positive control substance (4-OH ASDN), aromatase activity averaged 0.180±0.027 nmol·mg-protein<sup>-1</sup>·min<sup>-1</sup> at the lowest tested concentration (10<sup>-10</sup> M) and -0.008±0.004 nmol·mg-protein<sup>-1</sup>·min<sup>-1</sup> at the highest tested concentration (10<sup>-5</sup> M). The mean aromatase activity of the positive control (expressed as % full control activity) for each concentration tested across all three runs is presented in Table 7, along with the overall standard deviation and %CV. Inhibition response curves for the positive control from each run and the average of all runs are shown in Figure 1. These results were generally within the recommended ranges for the top of the curve, bottom curve, Hill slope, log IC<sub>50</sub>, and coefficient of variation for replicates of each concentration, with the following

b If the test compound was not soluble above  $10^{-6}$  M and the inhibition curve does not cross 50%, the chemical is typically determined to be un-testable in the aromatase assay.

exceptions: the overall %CV for the highest two concentrations (-45.2 and 415.9%) exceeded the acceptable limit of 15%, and the bottom of the curve in Runs 2 and 3 (-8.2 and -7.2, respectively) exceeded the acceptable range (-5 to +6).

TABLE 7. Effect of 4-OH ASDN and 2,4-D on Aromatase Activity (as percent of control) from Independent Runs <sup>a</sup>							
Chemical	Concen. (Log M)	# Runs	Overall Mean	Overall SD	Overall SEM	Overall %CV	
4-OH ASDN	-5	3	-4.49	2.03	1.17	-45.2	
(positive control)	-6	3	1.10	4.59	2.65	415.9	
	-6.5	3	16.86	2.40	1.39	14.2	
	-7	3	34.80	3.96	2.28	11.4	
	-7.5	3	61.98	5.07	2.93	8.2	
	-8	3	81.46	3.35	1.93	4.1	
	-9	3	92.09	1.65	0.95	1.8	
	-10	3	97.17	4.21	2.43	4.3	
2,4-D	-4	1	91.15	7.05	4.07	7.7	
	-4.5	3	92.20	7.34	4.24	8.0	
	-5	3	88.05	5.65	3.26	6.4	
	-6	3	91.77	4.17	2.41	4.5	
	-7	3	98.02	13.05	7.54	13.3	
	-8	3	105.55	9.09	5.25	8.6	
	-9	3	100.33	22.05	12.73	22.0	
	-10	3	104 10	11 19	6.46	10.7	


Values were calculated by the reviewers based on data provided on pages 66-68.

SD Standard Deviation

SEM Standard Error of the Mean

CV Coefficient of Variance

FIGURE 1. Inhibition Response Curves for 4-OH ASDN.

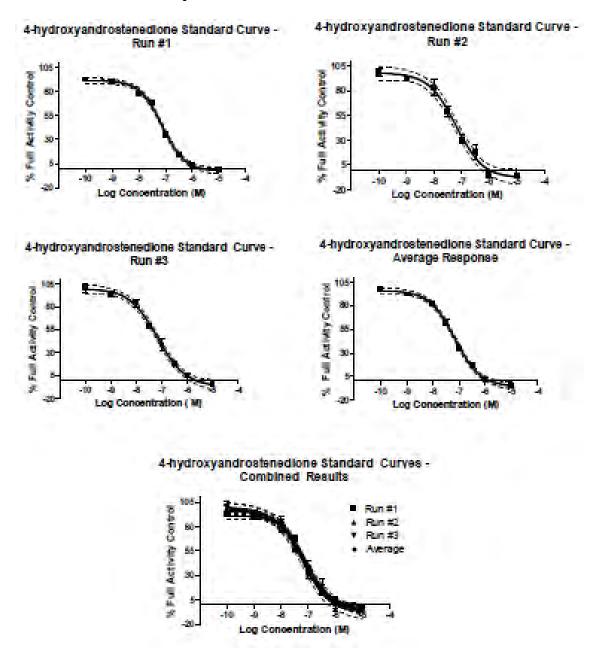

C. <u>TEST SUBSTANCE</u>: For 2,4-D, aromatase activity averaged 0.195±0.045 nmol·mg-protein<sup>-1</sup>·min<sup>-1</sup> at the lowest tested concentration (10<sup>-10</sup> M) and 0.168±0.025 nmol·mg-protein<sup>-1</sup>·min<sup>-1</sup> at the highest tested concentration (10<sup>-4</sup>). The mean aromatase activity of 2,4-D (expressed as %full control activity) for each concentration tested across all three runs is presented in Table 7 (above), along with the overall standard deviation and %CV. Inhibition response curves for 2,4-D from each run are shown in Figure 2, and the average inhibition response curve across all runs is shown in Figure 3. The overall %CV for the 10<sup>-9</sup> M concentration (22%) exceeded the acceptable limit of 15%.

FIGURE 2. Inhibition Response Curves for 2,4-D From Each Test Run.


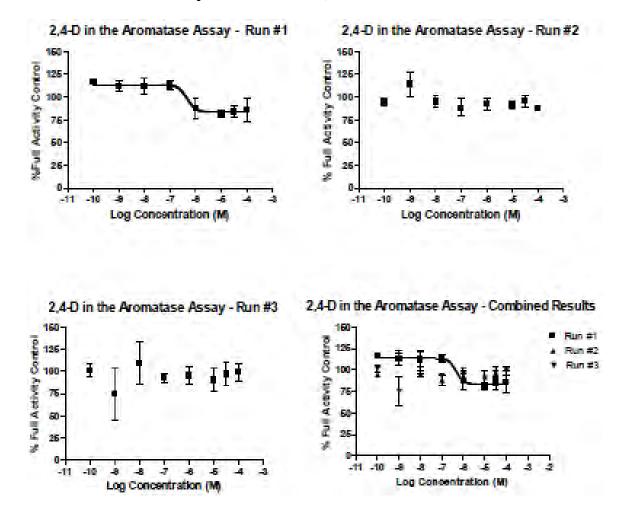


FIGURE 3. Mean Inhibition Response Curve for 2,4-D.

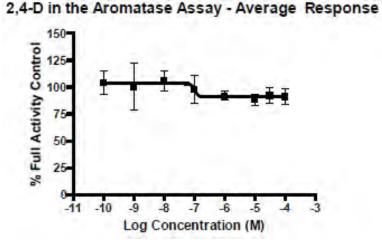

The data for 2,4-D were modeled for Run 1, and the goodness of fit (R<sup>2</sup>) value was 0.83; however, data from Runs 2 and 3 were not amenable to modeling. The average dose-

response curve indicated that the aromatase activity of the test material at concentrations ranging from  $10^{-10}$  M to  $10^{-4}$  M was essentially equivalent to the activity observed in the full activity controls.

For positive control curves, the best fit Hill slope and log IC<sub>50</sub> values across three independent runs were similar.

The effect of the positive control on inhibition of aromatase activity is presented in Table 8. For 4-OH ASDN, the estimated log IC<sub>50</sub> averaged -7.23 M and the slope was -0.92. Confidence in the mean log IC<sub>50</sub> for the positive control is high due to the small variation (<1% CV). Confidence in the mean slope is low due to the large variation (17.8% CV).

TABLE 8. Effect of 2,4-D on Aromatase Activity (as Percent of Control) From Independent Runs a						
Chemical	Run 1	Run 2	Run 3	Mean	SD	%CV
	Log IC <sub>50</sub> (M)					
2,4-D	NA	NA	NA	NA	NA	NA
4-OH ASDN	-7.2	-7.3	-7.2	-7.23	0.06	-0.80
		Sle	ope			
2,4-D	NA	NA	NA	NA	NA	NA
4-OH ASDN	-1.1	-0.88	-0.78	-0.92	0.16	-17.79

a Data were provided on page 32 of the study report. Mean, SD and %CV were calculated by the reviewers based on these data.

Based on the data from the average response curve and the criteria listed above in Table 6, the results support the conclusion that 2,4-D is a non-inhibitor of aromatase activity.

#### III. DISCUSSION AND CONCLUSIONS

- **A.** <u>INVESTIGATORS CONCLUSIONS</u>: The average response from the three independent runs with 2,4-D did not fit the four parameter regression model. Additionally, average aromatase activity for 2,4-D was similar to full activity controls at all concentrations tested. Therefore, 2,4-D is classified as a non-inhibitor of aromatase activity.
- **B.** <u>AGENCY COMMENTS</u>: Aromatase activity in the full activity controls ranged from 0.131 to 0.244 nmol·mg-protein<sup>-1</sup>·min<sup>-1</sup> for the three test runs, with a mean and standard deviation of 0.186±0.036 nmol·mg-protein<sup>-1</sup>·min<sup>-1</sup>. Activity in the background controls ranged from 7.31 to 11.51% of the full activity controls. The responses of the full activity controls were outside of the 90 to 110% range in the 2<sup>nd</sup> and 3<sup>rd</sup> replicates in the Run 1 (86.6 and 113.8%, respectively).

For the positive control substance (4-OH ASDN), aromatase activity averaged  $0.180\pm0.027$  nmol·mg-protein<sup>-1</sup>·min<sup>-1</sup> at the lowest tested concentration ( $10^{-10}$  M) and - $0.008\pm0.004$  nmol·mg-protein<sup>-1</sup>·min<sup>-1</sup> at the highest tested concentration ( $10^{-5}$  M). These results were within the recommended ranges for the top of the curve, bottom curve, Hill slope, log IC<sub>50</sub>, and coefficient of variation for replicates of each concentration, with the following

SD Standard Deviation

CV Coefficient of Variance

NA Not applicable. Values for 2,4-D were not suitable for modeling.

exceptions: the overall %CV for the highest two concentrations (-45.2 and 415.9%) exceeded the acceptable limit of 15%, and the bottom of the curve in Runs 2 and 3 (-8.2 and -7.2, respectively) exceeded the acceptable range (-5 to +6).

For 2,4-D, aromatase activity averaged  $0.195\pm0.045$  nmol·mg-protein<sup>-1</sup>·min<sup>-1</sup> at the lowest tested concentration  $(10^{-10}\,\text{M})$  and  $0.168\pm0.025$  nmol·mg-protein<sup>-1</sup>·min<sup>-1</sup> at the highest tested concentration  $(10^{-4})$ . The overall %CV for the  $10^{-9}\,\text{M}$  concentration (22%) exceeded the acceptable limit of 15%. The data for 2,4-D were modeled for Run 1, and the goodness of fit  $(R^2)$  value was 0.83. The average dose-response curve indicated that the aromatase activity of the test material at concentrations ranging from  $10^{-10}\,\text{M}$  to  $10^{-4}\,\text{M}$  was essentially equivalent to the activity observed in the full activity controls.

For 4-OH ASDN, the estimated log IC<sub>50</sub> averaged -7.23 M and the slope was -0.92. Confidence in the mean log IC<sub>50</sub> for the positive control is high due to the small variation (<1% CV), but confidence in the mean Hill slope is low due to the large variation (17.8% CV).

At  $10^{-4}$  M aromatase activity was 91% compared to the full activity controls. Since the average lowest portion of the activity response curve was greater than 75% activity, 2,4-D is classified as a non-inhibitor of aromatase activity up to the highest concentration tested ( $10^{-4}$  M).

#### C. STUDY DEFICIENCIES: None

# DATA EVALUATION RECORD

2,4-DICHLOROPHENOXY ACETIC ACID (2,4-D)

Study Type: OCSPP 890.1250, Estrogen Receptor Binding Assay

EPA Contract No. EP10H001452 Task Assignment No. 2-26-2012 (MRID 48614303)

Prepared for
Health Effects Division
Office of Pesticide Programs
U.S. Environmental Protection Agency
2777 South Crystal Drive
Arlington, VA 22202

Prepared by CSS-Dynamac Corporation 1910 Sedwick Road, Building 100, Suite B Durham, NC 27713

Primary Reviewer:	Signature:	wintell I She for
Michelle Sharpe-Kass, M.S.	Date:	3/15/2012
Secondary Reviewer:	Signature:	Lutenley
Scott D. Studenberg, Ph. D., D.A.B.T.	Date:	3/22/2012
Program Manager:	Signature:	Jack D. Eny
Jack D. Early, M.S.	Date:	3/27/2012
Quality Assurance:	Signature:	Jack D. Eury
<u>Jack D. Early, M.S.</u>	Date:	3/27//2012

This Data Evaluation Record may have been altered by the Health Effects Division subsequent to signing by CSS-Dynamac Corporation personnel.

The US EPA Endocrine Disruptor Screening Program (EDSP) Tier 1 screening battery is comprised of eleven screening assays intended to identify a chemical's likely endocrine bioactivity, i.e., its potential to interact with the estrogen, androgen, or thyroid (E, A, or T) pathways. The robustness of the Tier 1 battery is based on the strengths of each individual assay to identify potential endocrine bioactivity with complementary endpoints within the assay, where available, and redundancy across the battery. Thus, the results of each individual assay should not be considered in isolation but rather should be considered in the context of other assays in the battery as well as Other Scientifically Relevant Information (OSRI). In order to determine if a chemical has the potential to interact with the E, A or T pathways, a Weight of Evidence (WoE) evaluation of Tier 1 assay results, in combination with the findings in the OSRI, should be undertaken (refer to the WoE Document).

Primary Reviewer: Patience Browne, Ph.D.

Office of Science Coordination and Policy

Secondary Reviewer: <u>Greg Akerman, Ph.D.</u>

**Health Effects Division** 

Signature: Date:

Signature:

Date: 6/19/10

# DATA EVALUATION RECORD

**STUDY TYPE:** Estrogen Receptor Binding Assay Using Rat Uterine Cytosol (ER-RUC);

OCSPP 890.1250

<u>PC CODE</u>: 030001 <u>DP BARCODE</u>: D398640

TXR#: 0052104 CAS No.: 94-75-7

TEST MATERIAL (PURITY): 2,4-D (98.5% a.i.)

**SYNONYMS**: 2,4-dichlorophenoxy acetic acid

CITATION: LeBaron, M.J., Shisler, M.R., Visconti, N.R. (2011). Evaluation of 2,4-

Dichlorophenoxy Acetic Acid (2,4-D) in an *in vitro* Estrogen Receptor Binding Assay. Toxicology & Environmental Research and Consulting, Dow Chemical Company, Midland, MI. Laboratory Study ID.: 111121, October 27, 2011.

MRID 48614303. Unpublished.

SPONSOR: Industry Task Force II on 2,4-D Research Data, c/o McKenna Long & Aldridge

LLP, 1900 K Street NW, Washington, DC.

TEST ORDER #: CON-030001-1

EXECUTIVE SUMMARY: In an estrogen receptor (ER) binding assay (MRID 48614303) for 2,4-D (98.5%, Lot# 2006 2433 8006-USA), uterine cytosol from Sprague Dawley rats was used as the source of ER to conduct saturation binding and competitive binding experiments in this assay. The competitive binding experiment was conducted to measure the binding of a single concentration of [ $^3$ H]-17β-estradiol (1 nM) in the presence of logarithmic increasing concentrations of 2,4-D from  $10^{-11}$  to  $10^{-4}$  M, rather than  $10^{-10}$  to  $10^{-3}$  recommended in the test guideline. The justification for lowering the top test concentration was based on toxicokinetic data in the rat; concentrations higher than  $10^{-4}$  M were not relevant for testing in this assay as they are substantially above the inflection point for linear toxicokinetics. Ethanol was used as a solvent at a final concentration of <3%. The assay included 19-norethindrone as a weak positive control, octyltriethoxysilane as a negative control, and 17β-estradiol as the natural ligand reference material, and three independent runs were performed on separate days.

Summary data pertaining to the saturation binding experiment were reported separately in the study profile submitted by the test order recipient. The  $K_d$  for [ $^3H$ ]-17 $\beta$ -estradiol was 0.1032 nM and the  $B_{max}$  (nM) was 0.07097 for the prepared rat uterine cytosol used in these experiments. The  $K_d$  for the run was within the expected range of 0.03 to 1.5 nM, and the  $B_{max}$  was within the expected range of 10-150 fmol/100  $\mu$ g protein. The data produced a linear Scatchard plot.

In the competitive binding experiment, no precipitation was observed at any concentration tested. The mean specific binding in the presence of 2,4-D was >95% at 2,4-D concentrations of  $\leq 10^{-4}$  M in all three runs. The estimated mean log IC<sub>50</sub> and RBA was not calculated for 2,4-D as the percent binding inhibition did not reach 50% for any run.

The estimated mean log IC<sub>50</sub> for the natural ligand,  $17\beta$ -estradiol, and the weak positive control (19-norethindrone) was -9.0 and -5.5 M, respectively. The mean RBA was 0.034% for 19-norethindrone. All performance criteria were met for  $17\beta$ -estradiol, 19-norethindrone and octyltriethoxysilane.

2,4-D was tested over a concentration range that fully defined the top of the curve. The mean specific radioligand binding in the presence of 2,4-D was >95% at 2,4-D concentrations of  $\leq 10^{-4}$  M. Based on the results from the three runs, 2,4-D is classified as Not Interactive in the Estrogen Receptor Binding Assay.

The assay **satisfies** the EDSP Tier 1 Test Order requirements for an Estrogen Receptor Binding assay (OCSPP 890.1250).

**<u>COMPLIANCE</u>**: Signed and dated GLP, Data Confidentiality, and Quality Assurance statements were provided.

#### I. MATERIALS AND METHODS

## A. MATERIALS

1. Test Facility: Toxicology & Environmental Research and Consulting;

Dow Chemical

Location:Midland, MIStudy Director:M.R. Schisler

Other Personnel: M.J. LeBaron (Lead Scientist), N.R. Visconti (Research Biologist), B.B.

Gollapudi, (Technical Reviewer)

Study Period: September 14, 2011 - October 27, 2011

2. Test substance: 2,4-D

Description:Technical, Off-white PowderSource:NuFarm Americas, Inc.Lot/Batch #:2006 2433 8006-USA

**Purity:** 98.5%

**Solubility:** 315 mg/L in water, 30 mM in ethanol

Volatility:1.9 x 10-5 Pa at 25°CStability:2 year shelf lifeStorage conditions:AmbientCAS #:94-75-7

Molecular weight: 221.0

**Structure:** 

OH CI CI

**3.** Non-labeled ligand: 17β-estradiol

Supplier: Sigma, St. Louis MO

 Catalog #:
 E8875

 Batch #:
 098K1372

 Purity:
 100%

 CAS #:
 50-28-2

**1.** Radioactive ligand:  $[^3H]-17\beta$ -estradiol

Supplier: Perkin-Elmer, Boston MA

 Catalog #:
 NET517001MC

 Batch#:
 639068

Radiochemical purity: 597%
Specific activity: 162.9 Ci/mmol
Concentration of stock: 1.0 mCI/mL

5. <u>Positive control</u>: 19-norethindrone

Supplier: Sigma, St Louis, MO

 Catalog #:
 N4128

 Batch #:
 030M1359

 Purity:
 99%

 CAS #:
 68-22-4

6. Negative control: Octyltriethoxysilane

Supplier: Sigma, St Louis, MO

 Catalog #:
 440213

 Batch #:
 72596AMV

 Purity:
 98.58%

 CAS #:
 2943-75-1

7. Solvent/vehicle control: Ethanol

**Justification for choice of** None reported; 2,4-D is more soluble in water than ethanol

solvent:

Final Concentration: <3%

# B. METHODS

- 1. Preparation of Rat Uterine Cytosol (RUC): Trimmed uterine tissue from 85-100 day old female (Crl:CD(SD)) rats that were ovariectomized approximately 7-10 days prior to tissue harvest was purchased from Charles River (Wilmington, MA). Tissues were stored at −10° C until use (up to 6 months). The uteri were weighed, placed in ice-cold TEDG (Tris, EDTA, DTT, glycerol) + PMSF (phenylmethylsulfonyl fluoride) buffer and homogenized, followed by centrifugation for 10 min at 2500 × g at 4° C. Supernatant was transferred and centrifuged for 60 minutes at 105,000 × g, discarding the resulting pellets. Protein concentration of the cytosol was determined to be 3.573 mg/mL with the Pierce BCA method (Thermo Scientific Research Lab, Rockford, IL) using a protein kit compatible with DTT in the TEDG buffer. Cytosol was divided into aliquots (volume not reported) for immediate use or storage at −80° C for up to 90 days. The cytosol preparation was identified as Batch 2, prepared on 8/10/11.
- 2. <u>Saturation (radioligand) Binding Experiment</u>: A saturation binding experiment was conducted to demonstrate that the ER was present in adequate concentrations and had the appropriate affinity for the native ligand. A summary of the conditions for the saturation binding experiment are provided in Table 1.

TABLE 1. Summary of Conditions for Saturation Binding Experiment <sup>a</sup>					
Source of receptor		Rat uterine cytosol			
Concentration of radioligand	(as serial dilutions)	0.03-3.0 nM			
Concentration of non-labeled	ligand (100X [radioligand])	3.0-300 nM			
Concentration of receptor		Sufficient to bind 40.77% of radioligand at 0.03 nM <sup>b</sup>			
Temperature		~2-8 °C			
Incubation time		~16 hours			
Composition of assay buffer	Tris	10 mM (pH 7.4)			
	EDTA	1.5 mM			
Glycerol		10%			
Phenylmethylsulfonyl fluoride		1 mM			
	DTT	1 mM			

Data were not included in the study report, but were reported in the study profile submitted separately.

b This value was slightly higher than the suggested range in the guideline; however, all other values, including minimal ligand depletion, indicated acceptable performance in the assay.

The specific activity of the stock [³H]-17β-estradiol was not adjusted for decay over time on the day of the assay. Serial dilutions of radiolabeled estradiol in TEDG + PMSF buffer were prepared to achieve a final concentration of 0.03, 0.06, 0.08, 0.1, 0.3, 0.6, 1 and 3 nM. Solutions of non-labeled 17β-estradiol were prepared in a similar manner to achieve concentrations that were 100-fold greater than each respective radiolabeled concentration to result in final concentrations of 3, 6, 8, 10, 30, 60, 100 and 300 nM. For each batch of cytosol, the optimal protein concentration was determined by testing serial amounts of protein per tube, using 0.03 nM radiolabeled estradiol. The optimal protein concentration was determined to be 0.1191 mg protein/assay tube, which resulted in the binding of 40.77% of the total radioactivity added. This value was slightly higher than the suggested range in the guideline. Cytosolic protein used in this assay was thawed fresh for this experiment at ~4°C and maintained at ~4°C during the binding assay. Each run contained three concurrent replicates at each concentration, resulting in the 72 samples depicted in Table 2.

TABLE 2. Saturation Binding Experiment Run <sup>a</sup>				
Total binding <sup>b</sup>	Non-specific binding <sup>c</sup>	Radioligand alone <sup>d</sup>	Assay Components	
Tubes 1-24	Tubes 25-48	Tubes 49-72		
350 μL	300 μL		TEDG + PMSF buffer	
50 μL	50 μL	50 μL	[ <sup>3</sup> H]-17β-estradiol (8 serial dilutions) <sup>e</sup>	
	50 μL		Non-labeled 17β-estradiol (8 serial dilutions,	
			100x each respective labeled concentration) <sup>f</sup>	
100 μL	100 μL		Uterine cytosol (diluted to appropriate conc.)	
500 μL	500 μL	50 μL	Total volume in each assay tube	

- a Data were not included in the study report, but were reported in the study profile submitted separately.
- b Total binding =  $[^{3}H]$ -17β-estradiol bound to ER.
- c Non-specific binding =  $[^3H]$ -17 $\beta$ -estradiol and 100-fold greater non-labeled bound to ER.
- d Total [<sup>3</sup>H]-17β-estradiol alone for dpm determination at each concentration.
- e Final concentrations of  $[^{3}H]$ -17 $\beta$ -estradiol = 0.03, 0.06, 0.08, 0.1, 0.3, 0.6, 1, and 3 nM.
- f Final concentrations of non-labeled  $17\beta$ -estradiol = 3, 6, 8, 10, 30, 60, 100, and 300 nM.

Tubes were incubated with for ~16 hours at ~4°C. To separate bound from free estradiol, hydroxyapatite (HAP) slurry was added to each tube and vortexed (4 times with 5-minute intervals). Subsequently, the contents of each tube were washed three times as follows: TEDG +PMSF buffer was added, vortexed, centrifuged for 10 min at 1000 x g, and the supernatant decanted and discarded. After washing, ethanol was added to the HAP pellet remaining in each tube to extract the [3H]-17 $\beta$ -estradiol, followed by vortexing, and centrifugation for 10 min at 1000 x g. An aliquot of supernatant was radioassayed by scintillation counting. The temperature was maintained at approximately 4°C throughout the assay prior to extraction with ethanol. A total of two saturation binding runs on two batches of cytosol were performed with similar results. For the batch of cytosol used for the competitive binding assay, a single saturation run was performed.

3. <u>Competitive Binding Experiment</u>: A summary of the experimental conditions for the competitive binding experiment is presented in Table 3.

TABLE 3. Summary of Conditions for Competitive Binding Experiment <sup>a</sup>				
Source of receptor		Rat Uterine Cytosol		
Concentration of radioligand		1 nM		
Concentration of receptor		Sufficient to bind 6.75-7.15% of radioligand		
Concentration of test substance (as serial dilutions)		10 <sup>-11</sup> to 10 <sup>-4</sup> mM		
Temperature		4-8 °C		
Incubation time		16-20 hours		
Composition of assay buffer	Tris	10 mM (pH not reported)		
	EDTA	1.5 mM		
Glycerol		10% (v/v)		
Phenylmethylsulfonyl fluoride		1 mM		
	DTT	1 mM		

Data were obtained from pages 15, 21, 40, 42 and 44 of the study report.

The solubility of 2,4-D in ethanol was evaluated visually. On the day of the assay, the specific activity of the stock solution [³H]-17β-estradiol (not adjusted for decay over time) was diluted in TEDG + PMSF buffer to achieve a final concentration of 1 nM. For each batch of cytosol, the optimal protein concentration was determined by testing serial amounts of protein per tube, using 1.0 nM radiolabeled estradiol, until a concentration was reached that bound 6.75-7.15% of the total radioactivity added. Serial dilutions of the test substance, weak positive control (19-norethindrone), negative control (octyltriethoxysilane), and reference material (non-labeled 17β-estradiol) were prepared to achieve the concentrations shown in Table 4. Each assay consisted of three runs performed on separate days. Each run included three replicates of each test substance at each concentration, plus four blank and six samples of the master mix to determine full radioactivity, resulting in a total of 112 samples.

TABLE 4. Molar (M) concentrations in Competitive Binding Assay Run a, b					
	Positive control	Negative control	Reference Chemical		
2,4-D	19-norethindrone	Octyltriethoxysilane	Non-labeled 17β-estradiol		
Tubes 81-104 °	Tubes 33-56 °	Tubes 57-80 °	Tubes 11-32 and 105-112 °		
$10^{-10}$	$10^{-8.5}$	$10^{-10}$	Solvent control d		
$10^{-9}$	$10^{-7.5}$	$10^{-9}$	$10^{-11}$		
$10^{-8}$	$10^{-7}$	$10^{-8}$	$10^{-10}$		
$10^{-7}$	$10^{-6.5}$	$10^{-7}$	$10^{-9.5}$		
$10^{-6}$	$10^{-6}$	$10^{-6}$	$10^{-9}$		
$10^{-5}$	$10^{-5.5}$	$10^{-5}$	$10^{-8.5}$		
$10^{-4}$	$10^{-4.5}$	$10^{-4}$	$10^{-8}$		
	$10^{-4}$	$10^{-3}$	$10^{-7}$		

- a Data were obtained from pages 40-41 of the study report.
- b Each tube contains: 10μL of either the test substance, positive control, negative control, solvent control, or non-labeled 17β-estradiol; 390 μL of TEDG + PMSF buffer with [³H]-17β-estradiol; and 100 μL of uterine cytosol (with ER), for a total of 500 μL.
- c Each concentration of each chemical was run in triplicate, for a total of 96 tubes per run.
- d Solvent is ethanol

Tubes were incubated with gentle vortexing for 16-20 hours at  $4\pm2$  °C. To separate bound from free estradiol, hydroxyapatite (HAP) slurry was added to each tube and the tubes were vortexed (4 times with 5-minute intervals). Subsequently, the contents of each tube were washed three times as follows: TEDG+PMSF buffer was added, vortexed, centrifuged for 10 min at  $1000 \times g$ , and the supernatant decanted and discarded. Ethanol was then added to the HAP pellet remaining in each tube to extract the [ $^3$ H]-17 $\beta$ -estradiol, allowed to sit at room temperature for 15-20 min with vortexing (4 times with 5-minute intervals), and centrifugation for 10 min at  $1000 \times g$ . A portion of supernatant was radioassayed by scintillation counting. The temperature was maintained at  $4\pm2$  °C throughout the assay prior to extraction with ethanol.

C. <u>DATA ANALYSIS</u>: For the competitive binding experiment, total binding and non-specific binding data were modeled with a non-linear regression program [Graph Pad Prism v. 5.0 (GraphPad Software, Inc., La Jolla, CA)]. Nonlinear regression methods were used to fit a curve for 17β-estradiol, 19-norethindrone, octyltriethoxysilane, and 2,4-D data to the Hill equation with log IC<sub>50</sub> as a parameter to be estimated. Estimates of model parameters [e.g., log IC<sub>50</sub>, IC<sub>50</sub>, Hillslope, R<sup>2</sup>, and relative binding affinity (RBA)] were determined with Graph Pad Prism.

### 1. <u>Definitions</u>

**a.** <u>Classification of test material</u>: Classification of the test material is based on the average of three runs. Each run was first individually classified as follows:

Interactive = lowest point on the fitted curve within the range of the data is less than 50% (i.e., >50% of the radiolabeled estradiol has been displaced from the ER).

**Not interactive** = there are usable data points at or above  $10^{-6}$  M and either the lowest point on the fitted response curve within the range of the data is above 75% (i.e.,

<25% of the radiolabeled estradiol has been displaced from the ER) or a binding curve cannot be fitted and the lowest average percent binding among concentration groups in the data is above 75%.

Equivocal up to the limit of concentrations tested = there are no data points at or above a test chemical concentration of  $10^{-6}$  M and either a binding curve can be fit but  $\leq 50\%$  of the radiolabeled estradiol has been displaced from the ER or a binding curve cannot be fit and the lowest average percent binding among concentration groups in the data is  $\geq 50\%$ .

**Equivocal** = A run is classified as equivocal if it does not fall into any of the categories above

The categorical classification of each run was assigned a numerical value as follows:

Run Classification	Numerical Value
Interactive	2
Equivocal	1
Not interactive	0
Equivocal up to the limit of concentrations tested	"missing"

The values for each run were then averaged across runs and the chemical classified using the following ranges:

Test Material Classification	Numerical Range
Interactive	average ≥1.5
Equivocal	0.5≥ average <1.5
Not interactive	average < 0.5
Equivocal up to the limit of concentrations tested	"missing"

# b. <u>Descriptors for receptor binding:</u>

 $\mathbf{B}_{\text{max}}$ : maximum specific binding number (fmol ER/100 µg cytosol protein) measures the concentration of active receptor sites

K<sub>d</sub>: dissociation constant (nM), measures the affinity of the receptor for its natural ligand IC<sub>50</sub>: concentration of the test substance (M) at which 50% of the radioligand is displaced from the receptor

**Relative Binding Affinity (RBA %):** (IC<sub>50</sub> of 17 $\beta$ -estradiol  $\div$  IC<sub>50</sub> of test substance)  $\times$  100

#### II. RESULTS


A. SATURATION BINDING EXPERIMENT: Summary data pertaining to the saturation binding experiment were reported separately in the study profile submitted by the test order recipient. Individual data on saturation binding were not reported. Saturation binding experiment parameters are presented in Table 5. The K<sub>d</sub> for [<sup>3</sup>H]-17β-estradiol was 0.1032 nM and the B<sub>max</sub> (nM) was 0.07097 for the prepared rat uterine cytosol used in these experiments. The K<sub>d</sub> for the run was within the expected range of 0.03 to 1.5 nM. The B<sub>max</sub> was also within the expected range of 10-150 fmol/100 μg protein. Non-specific, specific

and total binding curves for [<sup>3</sup>H]-estradiol to the ER are shown in Figure 1. The data produced a linear Scatchard plot (Figure 2).

TABLE 5. Saturation Binding Experiment of 17β-estradiol with Estrogen Receptor from Rat Uterine						
Cytosol						
Parameter	Run 1	Run 2	Run 3	Runs 1-3		
R <sup>2</sup> (unweighted)	0.967	NR	NR	NR		
$B_{max}$ (nM)	0.07097	NR	NR	NR		
B <sub>max</sub> (fmol/100 μg protein)	59.28	NR	NR	NR		
$K_{d}$ (nM)	0.1032	NR	NR	NR		


a Data were not included in the study report, but were reported in the study profile submitted separately.


FIGURE 1. Binding of [ $^3$ H]-17 $\beta$ -Estradiol to the ER during the Saturation Binding Experiment.

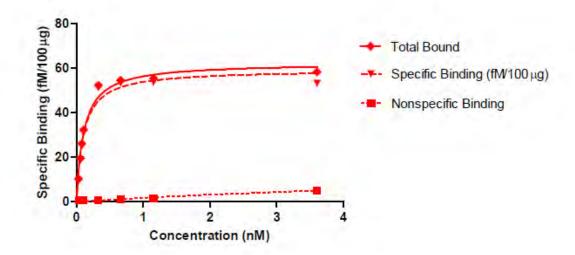

b Only a single run of the saturation binding experiment was conducted on the batch of cytosol used for this competitive binding experiment.

R<sup>2</sup> Goodness of fit for curve calculated for specific binding

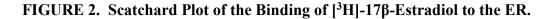

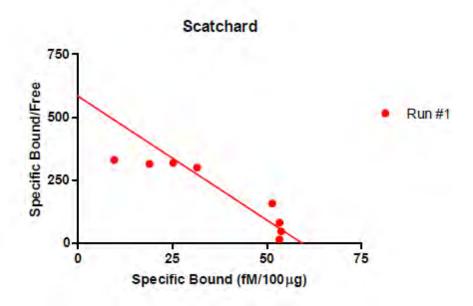

**B.** COMPETITIVE BINDING EXPERIMENT: The results from the three competitive binding experiments are summarized in Table 6 and presented graphically in Figures 3-5. No precipitation was observed at any concentration tested. The mean specific binding in the presence of 2,4-D was >95% at concentrations  $\leq 10^{-4}$  M in all three runs. The estimated mean log IC<sub>50</sub> and RBA was not calculated for 2,4-D as the percent binding inhibition did not reach 50% for any run.

The estimated mean log IC<sub>50</sub> for the natural ligand,  $17\beta$ -estradiol, and the weak positive control (19-norethindrone) was -9.0 and -5.5 M, respectively. The mean RBA was 0.034% for 19-norethindrone. As the lowest average percent binding in the presence of 2,4-D was  $\geq 95\%$  at concentrations up to  $10^{-4}$  M, 2,4-D is classified as not interactive (0) in this assay (Table 7).

TABLE 6. Co	TABLE 6. Competitive Binding Assay of 2,4-D with Estrogen Receptor from Rat Uterine Cytosol <sup>a</sup>					
Parameter		Run 1	Run 2	Run 3	$Mean \pm SE^b$	
r <sup>2</sup> (unweighted)	, 17β-estradiol	0.9990	0.9997	0.9995	0.9990-0.9997	
	19-norethindrone	0.9999	0.9997	0.9962	0.9962-0.9999	
	2,4-D	0.5875	NC	0.7841	0.5875-0.7841	
$Log IC_{50}(M),$	17β-estradiol	-8.948	-8.942	-8.961	$-8.951 \pm 0.006$	
	19-norethindrone	-5.462	-5.435	-5.534	$-5.477 \pm 0.030$	
	2,4-D	NC	NC	NC	NC	
$IC_{50}(M),$	17β-estradiol	$1.13 \times 10^{-9}$	$1.14 \times 10^{-9}$	$1.10 \times 10^{-9}$	$1.12 \times 10^{-9} \pm 0.01 \times 10^{-9}$	
	19-norethindrone	$3.45 \times 10^{-6}$	$3.76 \times 10^{-6}$	$2.92 \times 10^{-6}$	$3.33 \times 10^{-6} \pm 0.25 \times 10^{-6}$	
	2,4-D	NC	NC	NC	NC	
Log RBA (%),	19-norethindrone	-3.5	-3.5	-3.5	$-3.5 \pm 0.0$	
	2,4-D	NC	NC	NC	NC	
RBA (%),	19-norethindrone	0.033	0.031	0.032	$0.032 \pm 0.001$	
	2,4-D	NC	NC	NC	NC	

a Data were obtained from page 33 of the study report.


RBA (%) Relative binding affinity

NA Not applicable. r<sup>2</sup> is more appropriately expressed as a range, as opposed to a mean.

TABLE 7. Binding Classification of 2,4-D with Estrogen Receptor <sup>a</sup>					
Run	1	2	3	Mean <sup>c</sup>	Binding Classification d
Classification category value <sup>b</sup>	0	0	0	0	Not Interactive

a Data were obtained from pages 28-30 of the study report.

FIGURE 3. Percentage E2 Bound to the Estrogen Receptor in the Presence of Test Compound (Run 1).

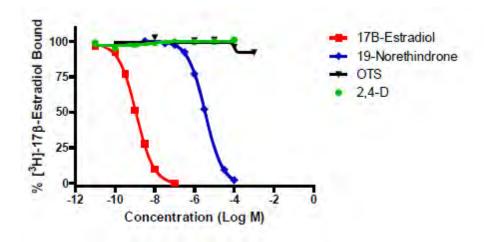

b The range is reported for  $r^2$ ; the mean  $\pm$  SE is reported for all other parameters.

r<sup>2</sup> Goodness of fit

b Classification category value: Interactive = 2; Equivocal = 1; Not interactive = 0; Equivocal up to the limit of concentrations tested ("missing", i.e., not included in calculation of mean).

c Mean of three runs expressed to the tenths place

d Interactive = mean  $\ge 1.5$ ; Equivocal =  $0.5 \le$  mean  $\le 1.5$ ; Not interactive = mean  $\le 0.5$ 

FIGURE 4. Percentage E2 Bound to the Estrogen Receptor in the Presence of Test Compound (Run 2).


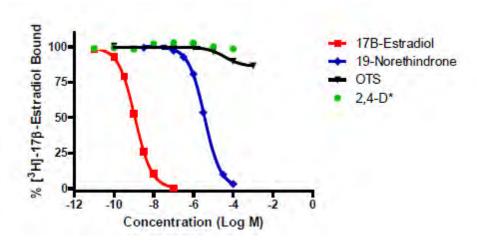


FIGURE 5. Percentage E2 Bound to the Estrogen Receptor in the Presence of Test Compound (Run 3).

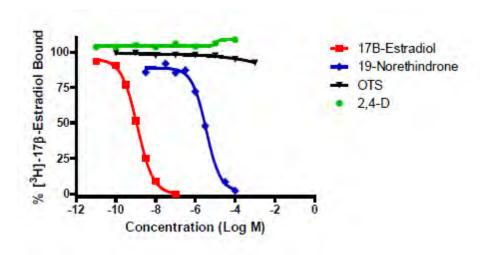

C. <u>PERFORMANCE CRITERIA</u>: To ensure that the competitive binding assay functioned properly, each run was evaluated using the following criteria (Table 8):

TABLE 8. Criterion <sup>a</sup>	Tolerance Limit(s)	Value	Yes	No
17β-estradiol fitted curve parameters				
Loge residual SD [Loge (SyX)]	≤2.35	-0.11 to 0.54	X	
Top (% binding)	94 to 111	95 to 99	X	
Bottom (% binding)	-4 to 1	−0.9 to −0.4	X	
Hill Slope $(\log_{10}(M)^{-1})$	−1.1 to −0.7	-1.08 to -0.98	X	
Weak Positive control (19-norethindrone) fitted curve paramete	ers <sup>b</sup>			
Log <sub>e</sub> residual SD	< 2.60	-0.52 to 1.10	X	
Top (% binding)	110 to 90	89 to 100		X
Bottom (% binding)	1 to -5	-1.2 to 1.7		X
Hill Slope $(\log_{10}(M)^{-1})$	-1.1 to -0.7	−1.2 to −0.99	X	
Solvent concentration				
Ethanol	≤3%	<3%	X	
Negative control (octyltriethoxysilane) does not displace more than 25% of [ ${}^{3}$ H]-17 $\beta$ -estradiol from the ER on average across all concentrations	≤25%	≤13.5	X	

Data were obtained from page 33 of the study report.

NA Not applicable

b The EPA Guideline does not define a set of tolerance limits for 19-norethindrone. Acceptance criteria were only defined for norethynodrel, which cannot be obtained commercially. The values reported were considered acceptable as they show 19-norethindrone to be an acceptable weak positive control.

Additionally, the curve for the reference material showed that increasing concentrations of unlabeled  $17\beta$ -estradiol displaced [ $^3$ H]- $17\beta$ -estradiol in a manner consistent with one-site binding, as indicated by a Hill slope of approximately -1.0. The recommended ranges of the weak positive control 19-norethindrone were occasional outside of the ranges for the top and bottom of curves established for norethynodrel. The fitted curve parameters demonstrate that 19-norethindrone is an acceptable weak positive control in the subject assay.

The percent binding of 2,4-D at this top plateau 98.8 - 103.5% was within 25 percentage points of the value for the lowest concentration of the estradiol standard 93.5 - 98.1%. Examination across the runs indicated consistency of the Hill slope, placement along the X-axis, and top and bottom plateaus.

### III. DISCUSSION AND CONCLUSIONS

- **A.** <u>INVESTIGATOR'S CONCLUSIONS</u>: Based on the combined responses in each of three independent estrogen receptor binding assays, it was determined that 2,4-D had no appreciable effect on the binding of the reference estrogen at any concentration, up to  $10^{-4}$ M. The results of the *in vitro* estrogen receptor binding assay using rat uterine cytosol indicate that, under the conditions of this study, 2,4-D was negative (not interactive) for estrogen receptor at concentrations up to  $10^{-4}$ M.
- **B.** AGENCY COMMENTS: The highest concentration of 2,4-D tested in this assay was 10<sup>-4</sup> M. The Test Guideline recommends testing up to 10<sup>-3</sup> M unless there is evidence of insolubility. The justification provided in the study report for lowering the top concentration was based on toxicokinetic data generated in the rat. The Test Guideline recommends testing up to 10<sup>-3</sup> M to adequately assess the potential for the test chemical to interact with the ER in mammalian as well as non-mammalian taxa. The justification for lowering the top concentration is inadequate since it does not apply to non-mammalian taxa.

Summary data pertaining to the saturation binding experiment were reported separately in the study profile submitted by the test order recipient; individual data on saturation binding were not reported. The  $K_d$  for  $[^3H]$ -17 $\beta$ -estradiol was 0.1032 nM and the  $B_{max}$  (nM) was 0.07097 for the prepared rat uterine cytosol used in these experiments. The  $K_d$  for the run was within the expected range of 0.03 to 1.5 nM, and the  $B_{max}$  was within the expected range of 10-150 fmol/100  $\mu$ g protein. The data produced a linear Scatchard plot (Figure 2).

The mean specific binding in the presence of 2,4-D was  $\geq$ 95% at concentrations  $\leq$ 10<sup>-4</sup> M in all three runs. The estimated mean log IC<sub>50</sub> and RBA was not calculated for 2,4-D as the percent binding inhibition did not reach 50% for any run

The estimated mean log IC $_{50}$  for the natural ligand, 17 $\beta$ -estradiol, and the weak positive control (19-norethindrone) was -9.0 and -5.5 M, respectively. The mean RBA was 0.034% for 19-norethindrone. All performance criteria were met for 17 $\beta$ -estradiol and octyltriethoxysilane. Examination of the data for the reference ligand and the weak positive control across the runs indicated consistency of the Hill slope, placement along the X-axis, and top and bottom plateaus.

- 2,4-D was tested over a concentration range that fully defined the top of the curve. The percent binding at this top plateau 98.8-103.5% was within 25 percentage points of the value for the lowest concentration of the estradiol standard 93.5-98.1%. Therefore, based on the combined responses in each of three independent estrogen receptor binding assay runs, it was determined that 2,4-D was not interactive with the estrogen receptor at concentrations up to  $10^{-4}$  M.
- **C. STUDY DEFICIENCIES:** The following deficiencies were noted that were not considered to have had an adverse impact on the results, interpretation or conclusions of this study:
  - The highest concentration of 2,4-D tested in this assay was 10<sup>-4</sup> M. The Test Guideline recommends testing up to 10<sup>-3</sup> M unless there is evidence of insolubility. The justification provided in the study report for lowering the top concentration was based on toxicokinetic data generated in the rat. The Test Guideline recommends testing up to 10<sup>-3</sup> M to adequately assess the potential for the test chemical to interact with the ER in mammalian as well as non-mammalian taxa.

### DATA EVALUATION RECORD

2,4-DICHLOROPHENOXYACETIC ACID (2,4-D)

Study Type: OCSPP 890.1300, Estrogen Receptor Transcriptional Activation

EPA Contract No. EP10H001452 Task Assignment No. 2-26-2012 (MRID 48614304)

Prepared for
Health Effects Division
Office of Pesticide Programs
U.S. Environmental Protection Agency
2777 South Crystal Drive
Arlington, VA 22202

Prepared by CSS-Dynamac Corporation 1910 Sedwick Road, Building 100, Suite B Durham, NC 27713

Primary Reviewer:	Signature:	milelle I She for
Michelle Sharpe-Kass, M.S.	Date:	2/19/2012
Secondary Reviewer:	Signature:	Lutenberg
Scott D. Studenberg, Ph.D., D.A.B.T.	Date:	3/26/2012
		Jack Q. Ewy
Program Manager:	Signature:	
Jack D. Early, M.S.	Date:	3/27/2012
		Jack Q. Ew4
Quality Assurance:	Signature:	
Jack D. Early, M.S.	Date:	3/27/2012

This Data Evaluation Record may have been altered by the Health Effects Division subsequent to signing by CSS-Dynamac Corporation personnel.

The US EPA Endocrine Disruptor Screening Program (EDSP) Tier 1 screening battery is comprised of eleven screening assays intended to identify a chemical's likely endocrine bioactivity, i.e., its potential to interact with the estrogen, androgen, or thyroid (E, A, or T) pathways. The robustness of the Tier 1 battery is based on the strengths of each individual assay to identify potential endocrine bioactivity with complementary endpoints within the assay, where available, and redundancy across the battery. Thus, the results of each individual assay should not be considered in isolation but rather should be considered in the context of other assays in the battery as well as Other Scientifically Relevant Information (OSRI). In order to determine if a chemical has the potential to interact with the E, A or T pathways, a Weight of Evidence (WoE) evaluation of Tier 1 assay results, in combination with the findings in the OSRI, should be undertaken (refer to the WoE Document).

Primary Reviewer: Ray Kent, Ph.D. Signature: Less On Ph.D.

Health Effects Division Date:

Secondary Reviewer: Minerva Mercado, Ph.D. Signature: Date: 6-30-15

Date: 6-30-/5 0
Template version 08/2011

### DATA EVALUATION RECORD

STUDY TYPE: Estrogen Receptor Transcriptional Activation (Human cell Line, HeLa-9903);

OCSPP 890.1300; OECD 455.

PC CODE: 030001 DP BARCODE: D398640

TXR#: 0052104 CAS No.: 94-75-7

**TEST MATERIAL (PURITY): 2,4-D (98.5% a.i.)** 

**SYNONYMS**: (2,4-Dichlorophenoxy) acetic acid

**CITATION:** LeBaron, M.J., Kah, H.L. (2011). Evaluation of 2,4-Dichlorophenoxy Acetic

Acid (2.4-D) in an *In Vitro* Estrogen Receptor Transcriptional Activation Assay in Human Cell Line HeLa-9903. Toxicology & Environmental Research and Consulting, Midland, MI. Laboratory Report No.: 111043, October 17, 2011.

MRID 48614304. Unpublished.

SPONSOR: Industry Task Force II on 2,4-D Research Data, c/o McKenna Long & Aldridge

LLP, 1900 K Street NW, Washington, D.C.

**TEST ORDER #: CON-03001-01** 

EXECUTIVE SUMMARY: In an estrogen receptor transcriptional activation assay (MRID 48614304), hERα-HeLa-9903 cells cultured *in vitro* were exposed to 2,4-D (98.5% a.i., Lot #2006 2433 8006 USA) at logarithmically increasing concentrations from 10<sup>-10</sup> to 10<sup>-4</sup> M in DMSO (0.1%) for 24 hours. A total of four separate runs were performed. Each run was performed using 96-well plates and each 2,4-D concentration was tested in triplicate (3 wells/plate). Cells were exposed to the test agent for approximately 24 hours to induce reporter (luciferase) gene products. Luciferase expression in response to activation of the estrogen receptor was measured upon addition of a luciferase substrate and detection with a luminometer with acceptable sensitivity.

2,4-D was tested up  $10^{-4}$  M based on solubility, cytotoxicity and *in vivo* toxicokinetic analysis. There were deviations from expected performance criteria for all of the four reference chemicals, but these deviations do not affect the interpretation of this study. The RPC<sub>Max</sub> was <0% for the first run, 8.8% for the second run, 5.3% for the third run, and 7.0% for the fourth run; the associated PC<sub>Max</sub> was  $10^{-4}$  M for runs 1-3 and  $10^{-5}$  M for run 4. Because the RPC<sub>Max</sub> < PC<sub>10</sub> in all assay runs, 2,4-D was considered negative for estrogen receptor transcriptional activation in this test system.

This assay satisfies the EDSP Tier 1 Test Order requirement for an Estrogen Receptor Transcriptional Activation assay (OCSPP 890.1300).

**COMPLIANCE:** Signed and dated GLP, Quality Assurance and Data Confidentiality statements were provided.

### I. MATERIALS AND METHODS

### A. MATERIALS

1. <u>Test Substance</u>: 2,4-D

Description: Technical, off-white powder
Source (Catalog #): Nufarm Americas (Not Reported)

**Lot/Batch #:** 2006 2433 8006 USA

**Purity:** 98.5%

**Solubility:** Up to 315 mg/L in water, soluble in DMSO up to 0.1M

Volatility: $1.9 \times 10^{-5}$  Pa at 25 °CStability:2 year shelf lifeStorage conditions:Ambient

Vapor pressure: 1.9 x 10<sup>-5</sup> Pa at 25°C

**CAS** #: 94-75-7

Structure:

OH OCI

### 2. Reference substances

17β-estradiol (strong estrogen; positive control)

Supplier: Sigma, St. Louis, MO

Catalogue and Batch #: Cat# E-8875, Lot# 079K0131

**Purity:** 100% **CAS #:** 50-28-2

 $17\alpha$ -estradiol (weak estrogen)

Supplier: Sigma, St. Louis, MO

**Catalogue and Batch #:** Cat # E-8750, Lot # 029K4116

**Purity:** ≥99.5% **CAS #:** 57-91-0

Corticosterone (negative compound)

Supplier: Sigma, St. Louis, MO

**Catalogue and Batch #:** Cat # C-2505, Lot # 010M2010

**Purity:** 100% **CAS #:** 50-22-6

 $17\alpha$ -methyltestosterone (very weak agonist)

Supplier: Sigma, St. Louis, MO

**Catalogue and Batch #:** Cat # M-7252, Lot # 060M1543V

Purity: 99% CAS #: 58-18-4

3. Vehicle(s)

Solvent: DMSO, Sigma-Aldrich, Cat # 276855

Solvent control 0.1%

(final concentration):

### B. METHODS

- 1. <u>Cell Culture</u>: Stably-transfected hERα-HeLa-9903 cells were obtained from the Japanese Collection of Research Bioresources Cell Bank and were verified to be free of mycoplasma infection by ATCC (method not reported). Cells were maintained in Eagles Minimum Essential Medium without phenol red, supplemented with 60 mg/L kanamycin and 10% dextran-coated charcoal-treated fetal bovine serum (DCC-FBS; Hyclone Laboratories, Inc., Logan, Utah Lot# not reported), in an incubator under 5% CO<sub>2</sub> at 37° C. Upon reaching 75-90% confluence, cells were subcultured at least twice prior to exposure to the test material.
- 2. Transcriptional Activation Assays: For each test, cells were plated at a density of 1 1.5×10<sup>4</sup> cells/100 μL medium/well in a 96-well plate and allowed to attach for at least 3 hours. Growth media was replaced with media containing serial log dilutions of 2,4-D in DMSO (0.1% total final concentration). Cells were incubated for approximately 24 hours at approximately 37° C. Cytotoxicity was determined by a modified MTT cell viability assay, when stock solutions were diluted with treatment media, after addition to the cell culture plate and after the 24-hour treatment period. Transcriptional activation of the estrogen receptor (firefly luciferase activity) was determined using a standard assay kit (Promega, Madison, WI). Chemiluminescence was measured using a Packard TopCount NXT luminescence counter with a detection limit of ~5 cells in 100μL medium.
  - a. <u>Preliminary Test</u>: A preliminary test evaluating concentrations ranging from  $10^{-4}$  to  $10^{-10}$  M was conducted to determine the appropriate concentration range and to determine concentrations resulting in insolubility and/or cytotoxicity.
  - **Proficiency Chemicals:** It was stated that laboratory validation assays with 10 proficiency chemicals were performed to confirm the responsiveness of the ER transcriptional activation assay. These non-GLP unpublished results were reported to have demonstrated laboratory proficiency.
  - **c.** <u>Reference Chemicals</u>: To ensure the stability of the response from the cell line, six concentrations of each of the following reference chemicals were included on each plate in the current assay, along with the test chemical:

Reference Chemical	CAS No.	Concentration Range	Class
17β-estradiol	50-28-2	$10^{-14}$ to $10^{-8}$	Strong estrogen
17α-estradiol	57-91-0	$10^{-12}$ to $10^{-6}$	Weak estrogen
Corticosterone	50-22-6	$10^{-10}$ to $10^{-4}$	Negative compound
17α-methyltestosterone	58-18-4	$10^{-11}$ to $10^{-5}$	Very weak agonist

3. <u>Data analysis</u>: To obtain the relative transcriptional activity to the 1 nM E2 positive control (PC), the luminescence signals from the concurrent plate were analyzed by subtracting the mean value of the vehicle control from each well value to normalize the data; each normalized value was then divided by the mean value of the normalized PC. The resulting value was multiplied by 100 in order to express relative transcriptional activity as a percentage of the PC. Graph Pad Prism v. 5.0 (GraphPad Software, Inc., La Jolla, CA) was used to calculate the EC<sub>50</sub>, PC<sub>10</sub>, PC<sub>50</sub>, RPC<sub>Max</sub>, and PC<sub>Max</sub> for 2,4-D when applicable. The test material was defined as negative for inducing estrogen receptor transcriptional activation if the RPC<sub>Max</sub> < PC<sub>10</sub> in at least 2 of 2 (or 2 of 3) runs. Log EC<sub>50</sub> and Hill slope

values are calculated only if a positive response is observed. Coefficients of variation (CV) were calculated for the luminescence data triplicates. Concentrations showing >20% cytotoxicity or evidence of insolubility were excluded from analyses.

### 4. **Definitions**

- $EC_{50}$  = concentration of agonist that induces a response halfway between the baseline (bottom) and maximum (top) response
- $PC_{10}$  = concentration of a test chemical at which the response is 10% of the response induced by the positive control (E2 at 1 nM) in each plate
- $PC_{50}$  = concentration of a test chemical at which the response is 50% of the response induced by the positive control (E2 at 1 nM) in each plate
- $RPC_{Max}$  = maximum level of response induced by a test chemical, expressed as a percentage of the response induced by the positive control (1 nM E2) on the same plate

 $PC_{Max}$  = concentration of a test chemical inducing the  $RPC_{Max}$ 

### II. RESULTS

**A.** PRELIMINARY TEST: 2,4-D was relatively non-toxic and freely soluble; therefore the highest concentration was set at 10<sup>-4</sup> M based on the *in vivo* 2,4-D plasma concentrations that exceeded the threshold for linear toxicokinetics. These data were provided in an appendix to the report. No solubility or cytotoxicity issues were noted at the concentrations tested. Based on these results, logarithmically increasing concentrations from 10<sup>-10</sup> to 10<sup>-4</sup> M were selected for the assay.

TABLE 1. Preliminary	Test for Solubility, Cyto	toxicity, and Concentration-Selection for 2,4-Da
Concentration (M)	% Viability	Comments
$10^{-4}$	119.8	Toxicokinetic-derived limit dose
$10^{-5}$	123.6	
$10^{-6}$	133.1	
$10^{-7}$	132.0	
$10^{-8}$	115.1	
$10^{-9}$	130.4	
$10^{-10}$	117.6	
E2 at 1 nM	130.8	
$VC^b$	100.0	

Data were obtained from page 51 of the study report.

### **B.** Positive and Negative Reference Chemicals

1. <u>Proficiency Chemicals</u>: The responsiveness of cells to the required proficiency chemicals was not reported, but it was stated that the proficiency validation assays were conducted, and proficiency was demonstrated.

b Vehicle Control

TABLE 2. Proficiency Chemicals <sup>a</sup>		
Compound	<b>Expected Response</b>	Lab Response
Diethylstilbestrol	Positive	Not reported
17α-Ethynyl estradiol	Positive	Not reported
Hexestrol	Positive	Not reported
Genistein	Positive	Not reported
Estrone	Positive	Not reported
Butyl paraben	Positive	Not reported
1, 3, 5-Tris(4-hydroxyphenyl)benzene	Positive	Not reported
Dibutyl phthalate	Negative	Not reported
Atrazine	Negative	Not reported
Corticosterone	Negative	Not reported

2. Reference Chemicals: Values derived from the concentration response curve (e.g., log PC<sub>50</sub>, log PC<sub>10</sub>, log EC<sub>50</sub>, Hill slope) for the four concurrently run reference materials are included in Table 3. There were deviations from expected performance criteria for all of the four reference chemicals. The log EC<sub>50</sub> was higher than the expected value and the Hill slope lower than the expected value for 17β-estradiol. The log PC<sub>50</sub> and log PC<sub>10</sub> for 17α-estradiol and 17α-methyltestosterone were lower than expected, as was the Hill slope for 17α-estradiol. For 17α-methyltestosterone, the mean RPC<sub>Max</sub> was 90.8% for the first run, 152.8% for the second run, 83.3% for the third run, and 67.4% for the fourth run. For corticosterone, the mean RPC<sub>Max</sub> was 6.3% for the first run, 35.5% for the second run, 7.0% for the third run, and 5.4% for the fourth run. The deviations from the expected values do not negatively affect the interpretation of this study. Performance criteria values below the validated ranges usually indicate increased sensitivity of the assay compared to validation experiments.

TABLE 3. Performan	nce Criteria for Refe	erence Cho	emicals <sup>a</sup>				
Reference Chemical	A ( . l. l D		Val	lues		Acce	ptable
Parameter	Acceptable Range	Run 1	Run 2	Run 3	Run 4	Yes	No
17β-estradiol							
Log PC <sub>50</sub>	-11.4 to -10.1	-10.8	-11.1	-10.9	-10.7	X	
Log PC <sub>10</sub>	<-11	-11.9	-12.6	-12.7	-12.5	X	
Log EC <sub>50</sub>	-11.3 to -10.1	-9.4	-10.4	-10.7	-10.6		Run 1
Hill Slope	0.7 to 1.5	0.3	0.3	0.7	0.8	X	
Test range	$10^{-14}$ to $10^{-8}$ M		$10^{-14}$ to	10 <sup>-8</sup> M		X	
17α-estradiol							
Log PC <sub>50</sub>	−9.6 to −8.1	-9.3	-10.1	-9.0	-9.0	X	
Log PC <sub>10</sub>	−10.7 to −9.3	-10.6	-10.9	-10.6	-10.9	X	
Log EC <sub>50</sub>	−9.6 to −8.4	-8.7	-9.6	-8.8	-8.6	X	
Hill Slope	0.9 to 2.0	0.5	1.1	0.7	0.6	X	
Test range	10 <sup>-12</sup> to 10 <sup>-6</sup> M		$10^{-12}$ to	10 <sup>-6</sup> M		X	
Corticosterone							
Test range	10 <sup>-10</sup> to 10 <sup>-4</sup> M		$10^{-10}$ to	10 <sup>-4</sup> M		X	
17α-methyltestosterone							
Log PC <sub>50</sub>	−6.0 to −5.1	-7.0	-8.8	-7.5	-6.5	X	
Log PC <sub>10</sub>	-8.0 to -6.2	-8.8	-9.8	-9.2	-9.2	X	
Test range	10 <sup>-11</sup> to 10 <sup>-5</sup> M		$10^{-11}$ to	10 <sup>-5</sup> M		X	

<sup>&</sup>lt;sup>a</sup> Data were obtained from page 27 of the study report.

### C. <u>DEFINITIVE ASSAY</u>

1. <u>Vehicle and Positive Controls</u>: Data for the vehicle and positive controls are included in Table 4. The overall mean TA value for the vehicle control was 1513-3035 arbitrary light

units, and the overall mean TA value for the positive control was 12128-19639 (not reported for Run 4). The induction for the positive control ranged from 4.7- to 8.2-fold. The mean normalized value for the positive control was 10912-16604. The PC<sub>50</sub> (50% of the maximum response) for E2 in this assay is 6212.5-9819.5 and the PC<sub>10</sub> (10% of the maximum response) is 1243-1964.

TABLE	4. Transcr	iptional Ac	tivation (TA)	Response of	Vehicle and Posit	ive Control <sup>a</sup>	
Sample	Vehicle	Control		Positive Con	trol <sup>b</sup>	Normalized Posi	tive Control b
Runs	Mean	SD	Mean	SD	Fold Induction <sup>c</sup>	Mean	SD
1	2743	248	15564	1604	5.7	12821	NR
2	1513	290	12425	1958	8.2	10912	NR
3	3035	271	19639	1116	6.5	16604	NR
4	2603	NR	12128 <sup>d</sup>	915 <sup>d</sup>	4.7 <sup>d</sup>	12158	NR

- Data were obtained from page 47, 49 and 53 of the study report.
- b Positive control was 17β-estradiol (E2) at 1 nM.
- <sup>c</sup> Fold-induction = (mean TA of PC)/(mean TA of VC)
- d Calculated by the reviewer from data on page 49 of the study report
- NR Not reported
- 2. <u>Test Material</u>: Relative (to the PC) transcriptional activation at each concentration of the test chemical during the three assay runs is presented in Table 5. The concentration-response bar graphs depicting fold induction of relative transcriptional activation are presented in Figure 1 below. The RPC<sub>Max</sub> was <0% for the first run, 8.8% for the second run 5.3% for the third run and 7.0% for the fourth run; the associated PC<sub>Max</sub> was  $10^{-4}$  M for Runs 1-3 and  $10^{-5}$  M for Run 4. Because the RPC<sub>Max</sub> < PC<sub>10</sub> in all four runs, 2,4-D was considered negative for estrogen receptor transcriptional activation in this test system.


TABLE 5. Relat	tive Transcri	iptional A	ctivation (R	RTA) of 2,4	1-D <sup>a</sup>			
Parameter			RTA (mean	n ± SD); %	of Positive C	ontrol (PC)		
	Ru	n 1	Ru	n 2	Ru	n 3	Ru	n 4
Conc. (M)	Mean	SD	Mean	SD	Mean	SD	Mean	SD
$10^{-4}$	-7.2	3.8	8.8	1.0	5.3	1.1	4.4	4.8
$10^{-5}$	-4.6	4.1	3.3	4.1	1.7	1.0	7.0	3.6
10-6	-10.1	0.9	-0.6	1.8	-6.6	0.4	2.4	0.3
$10^{-7}$	-8.9	3.4	-2.1	2.0	-6.0	3.6	-1.4	2.3
10-8	-8.8	5.0	-2.5	1.4	-5.8	5.4	-1.1	2.7
10-9	-6.2	4.8	-3.2	0.9	-7.5	4.7	0.1	1.9
$10^{-10}$	-1.9	3.6	-0.8	1.2	-2.5	3.8	4.3	1.1
Log EC <sub>50</sub> <sup>b</sup>	N.	A	N	ſΑ	N	A	N	A
Hill Slope <sup>b</sup>	N.	A	N	A	N	A	N	A
RPC <sub>Max</sub>	0	b	8	.8	5	.3	7.	.0
PC <sub>Max</sub>	N.	A	10	)-4	10	)-4	10	)-5
PC50	N.	A	N	A	N	A	N	A
PC <sub>10</sub>	N.	A	N	A	N	A	N	A

Data were obtained from page 54 of the study report.

NA Not Applicable

b Value of RPCMax less than or equal to 0

Figure 1. Fold Induction of Relative Transcription Activation (RTA) of 2,4-D Compared to the Positive Control.

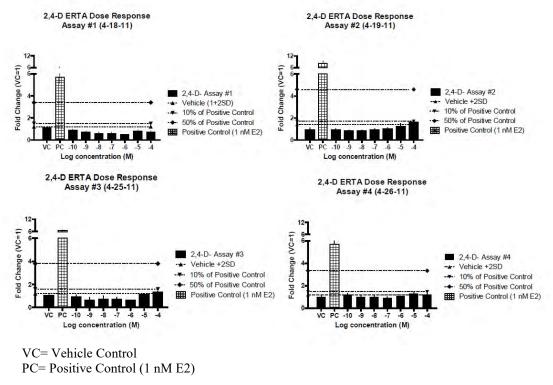

3. Performance Criteria: The results of the laboratory proficiency test were not reported. There were minor deviations from expected performance criteria for all of the four reference chemicals. The Log EC<sub>50</sub> was higher than the expected value and the Hill slope lower than the expected value for 17β-estradiol. The log PC<sub>50</sub> and log PC<sub>10</sub> for 17α-estradiol and 17α-methyltestosterone were lower than expected, as was the Hill slope for 17α-estradiol. The minor deviations from the expected values do not negatively affect the interpretation of this study. Mean luciferase activity was 4.7 to 8.2-fold. The fold-induction corresponding to the PC<sub>10</sub> of the concurrent PC was greater than 1+2 SDs of the vehicle control on all plates. Variability was minimal, and results were reproducible among runs.

### III. DISCUSSION AND CONCLUSIONS

- A. <u>INVESTIGATORS' CONCLUSIONS</u>: Based on the combined responses in each of four independent estrogen receptor transactivation assays, it was determined that 2,4-D treatment did not result in ER-mediated transcriptional activation at any concentration, including the toxicokinetic-derived limit, in this assay system (10<sup>-4</sup> M). The highest concentration of 2,4 D was based on *in* vivo toxicokinetic analyses; concentrations higher than 10<sup>-4</sup> M were not considered relevant for use in this assay as they are substantially above the inflection point for linear toxicokinetics.
- **B.** <u>AGENCY COMMENTS</u>: 2,4-D was tested up to the toxicokinetic-derived limit, 10<sup>-4</sup>M. The laboratory proficiency assays were not reported. The only deviation was a weak positive response in one run in response to corticosterone. There were minor deviations from expected performance criteria for all of the four reference chemicals, but these minor

deviations do not negatively affect the interpretation of this study. The RPC<sub>Max</sub> was <0% for the first run, 8.8% for the second run 5.3% for the third run and 7.0% for the fourth run; the associated PC<sub>Max</sub> was  $10^{-4}$  M for Runs 1-3 and  $10^{-5}$  M for Run 4. Because the RPC<sub>Max</sub> < PC<sub>10</sub> in all four runs, 2,4-D was considered negative for estrogen receptor transcriptional activation in this test system.

### **C. STUDY DEFICIENCIES:** The following deficiencies were noted:

- The laboratory proficiency assays were not included in the study report, but were sent in an accompanying study profile template.
- 2,4-D was not tested up to the limit dose, but rationale was provided to support the concentrations tested.

Data Requirement: EPA DP Barcode 388580

OECD Data Point

**EPA MRID** 48317001 EPA Guideline 890.1350

Fish Short-Term Reproduction Assay

Test Material: 2,4-Dichlorophenoxyacetic acid Purity (%): 98.6%

Common Name 2,4-D **IUPAC** Chemical Name

CAS Name

CAS No. 94-75-7 Synonyms 2,4-D acid EPA PC Code 030001

Digitally signed by PATIENCE BROWNE

DN: c=US, o=U.S. Government, ou=USEPA, ou=Staff, cn=PATIENCE BROWNE, dnQualifier=0000048202

Digitally signed by ROBIN STERNBERG

STERNBERG
DN: c=US, o=U.S. Government,
ou=USEPA, ou=Staff, cn=ROBIN
STERNBERG,
dnQualifier=0000039126
Date: 2015.06.01 13:17:23 -04'00'

Date: 2015.06.03 15:14:20 -04'00'

Primary Reviewer: Patience Browne Signature:

USEPA/OCSPP/OSCP Date: 06/25/2012

Signature: No longer with EPA Additional Reviewer: Alicia Korol

USEPA/OCSPP/OPP/EFED/ERB1 Date: 06/14/2011

**AMY** Additional Reviewer: Amy Blankinship Signature: **BLANKINSHIP** 

USEPA/OCSPP/OPP/EFED/ERB3 Date: 12/28/2012

Final Additional Reviewer: Robin Sternberg Signature:

USEPA/OCSPP/OPP/EFED/ERB1 Date: 05/27/2015

Date Evaluation Completed: 05/27/2015

CITATION: Marino, T.A., K.K. Coady, L.K. Sosinski, J. Thomas. 2010, DICHLOROPHENOXYACETIC ACID: A FISH SHORT-TERM REPRODUCTION ASSAY USING THE FATHEAD MINNOW, *Pimephales promelas*, Toxicology and Environmental Research and Consulting, The Dow Chemical Company. and Midland, Michigan 48674, Laboratory ID: 101026, Industry Task Force II on 2,4-D Research Data c/o McKenna Long & Aldridge LLP, Washington, D.C, Completed on 06 December 2010.

The US EPA Endocrine Disruptor Screening Program (EDSP) Tier 1 screening battery is comprised of eleven screening assays intended to identify a chemical's likely endocrine bioactivity, i.e., its potential to interact with the estrogen, androgen, or thyroid (E, A, or T) pathways. The robustness of the Tier 1 battery is based on the strengths of each individual assay to identify potential endocrine bioactivity with complementary endpoints within the assay, where available, and redundancy across the battery. Thus, the results of each individual assay should not be considered in isolation but rather should be considered in the context of other assays in the battery as well as Other Scientifically Relevant Information (OSRI). In order to determine if a chemical has the potential to interact with the E, A or T pathways, a Weight of Evidence (WoE) evaluation of Tier 1 assay results, in combination with the findings in the OSRI, should be undertaken (refer to the WoE Document).

Disclaimer: The guideline recommendations in this DER template are offered as a general reference to aid in preparation of the DER. The purpose of these recommendations is not to serve as substitute for the Test Guidelines, nor to provide any guidance on how the study should be conducted.

**EXECUTIVE SUMMARY** 

The 21-day short-term reproduction assay of 2,4-D with fathead minnow (Pimephales promelas) was conducted

under flow-through conditions. Adult fish (20 spawning groups; 2 males and 4 females in each group; 4

groups/treatment; ca. 6 months old) were exposed to 2,4-D (98.6% purity) at nominal concentrations of 0

(negative control), 0.400, 4.00, 40.0, and 100 mg a.i./L concentrations with corresponding mean-measured

concentrations of <0.10 (<LOQ, negative control), 0.245, 3.14, 34.0, and 96.5 mg a.i./L, respectively. The

test system was maintained at 24.5 to 25.2°C and a pH of 7.02 to 7.76.

The single mortality observed during the assay occurred in a female in the high treatment group. There were

no significant differences for male or female body weight or length relative the negative control. At test

termination, observations of secondary sex characteristics were observed in the negative control and treated

groups; no treatment-related effects were reported. Clinical signs included the loss of an eye, ascites, and

scoliosis (bent tail) which were observed in single fish in the negative control or treatment groups.

Spawning occurred in the negative control at least every 4 days in 3 of the 4 replicates, and mean fecundity

was 29.8 eggs/female/day/replicate; fertility in the negative control was 96.9%. Fecundity was significantly

decreased (Jonckheere-Terpstra; p<0.05) by 34% in the high treatment group with a non-significant (p>0.05)

concentration-dependent trend of decreased fecundity in the lower treatment groups compared to the negative

control. There were no significant differences for fertility between the 2,4-D treatments and the negative control.

There were no significant differences (p>0.05) between the 2,4-D treatment groups and the negative control

for male or female gonado-somatic index (GSI) or plasma vitellogenin (VTG). There was also no significant

difference (p>0.05) for male nuptial tubercle scores; no tubercles were observed for females. No apparent

treatment-related histopathological effects were observed in males and females. Although not concentration-

dependent, an increase in the number of female ovaries that were observed as Stage 2 (compared to Stage

3 or 4) was reported for the 2,4-D treatments compared to the negative control. Plasma sex steroid

concentrations were not reported.

The performance and validity criteria were met in this study with the exception that the coefficient of variation

(CV) for the mean-measured concentration of the lowest treatment group was 43%, exceeding the guideline

criterion of <20%. This deviation did not impact the interpretation of the study.

This assay satisfies the EDSP Tier 1 Test Order requirements for a Fish Short-Term Reproduction Assay (OCSPP Guideline 890.1350).

### Results Synopsis:

Test organism age at test initiation: 6 months

Mean body weight at test initiation: 3.7 g for males, 1.9 g for females

Mean length at test initiation: Not reported

Test type: Flow-through

## Data Evaluation Record on the Fish Short-Term Reproduction Assay with 2,4-D

### EPA MRID Number 48317001

Table 1: Summary of Reproductive and HPG Effects<sup>1,2</sup> in the Fish Short-Term Reproduction Assay (FSTRA) with 2,4-D.

Treatment			Tubercle Score	Score	5	GSI	Gonada	Gonadal Histo.	Plasma VTG	VTG	Plasma T	па Т	Plasma E2	ia E2
(mg a.i./L)	Focupality	Fert.												
[mean-	contain	Success	Σ	ш	Σ	ш	Σ	ш	Σ	ш	Σ	ш	Σ	ш
measured]														
0.245	No	No	No	No	No	No	oN	No	No	No	NA	NA	NA	NA
3.14	No	No	No	%	No	No	oN	No	No	No	NA	AN	NA	Ą
34.0	No	No	No	%	N <sub>o</sub>	No	oN	No	No	No	NA	AN	NA	Ą
96.5	Yes	No	No	No	No	No	No	No	No	No	NA	AN	NA	NA

Abbreviations: <sup>Conc.</sup> Concentration. <sup>Diff.</sup> Difference. <sup>E2</sup> 17β-estradiol. <sup>F</sup> Female. <sup>Fert.</sup> Fertilization. <sup>GSI</sup> Gonado-Somatic Index. <sup>Histo</sup> Histopathology.

 $^{\rm M}$  Male.  $^{\rm NA}$  Not applicable.  $^{\rm T}$  Testosterone.  $^{\rm VTG}$  Vitellogenin.

A "yes" indicates a significant difference based on comparison to the negative (clean water) control, unless otherwise specified.

<sup>2</sup> The criteria for significance are described in the Reviewer's Analysis and Statistical Verification sections of the DER. Conclusions regarding histopathology may be heavily weighted by the expert opinion of a board-certified pathologist.

Page 5 of 86

### I. MATERIALS AND METHODS

Guideline Followed:

This study was conducted in accordance with the Endocrine Disruptor Screening Program Guidelines, OCSPP (form. OPPTS) 890.1350: Fish Short-Term Reproduction Assay, Environmental Protection Agency (EPA) 740-C-09-007, October 2009. Supporting guidance documents used include: Organization for Economic Cooperation and Development (OECD) Guideline for Testing of Chemicals No. 229: "Fish Short Term Reproduction Assay", 2009; "A Short-term Test Method for Assessing the Reproductive Toxicity of Endocrine-Disrupting Chemicals Using the Fathead Minnow (Pimephales promelas)", EPA/600/R-01/067, 2002; and "Guidance Document on Aquatic Toxicity Testing of difficult Substances and Mixtures", OECD Series on Testing and Assessment. No. 23. 2000. Deviations from the OCSPP 890.1350 include the following:

- The coefficient of variation (CV) for the mean-measured concentration exceeded 20% in the 0.245 mg a.i/L treatment group.
- 2. Ammonia levels were detected at 0.24 mg/L, exceeding guideline recommendations of ≤1 μg/L.
- 3. Organic carbon levels were not reported in the routine laboratory dilution water solute analysis.
- 4. OCSPP guidelines recommend chlorine levels at <10  $\mu$ g/L; chlorine in the dilution water was not detectable at a level of 0.20 mg/L and it is unknown whether levels exceeded 10  $\mu$ g/L.

These deviations do not impact the interpretation of the study.

Compliance:

Signed and dated GLP, Quality Assurance, and No Data Confidentiality statements were provided. All phases of this study were conducted in accordance with the GLP Principles of the USEPA – FIFRA GLPs Title 40 CFR, Part 160 – Federal Insecticide Fungicide Rodenticide Act (FIFRA), Good Laboratory Practice Standards, Final Rule; OECD Series on principles of Good Laboratory Practice and Compliance Monitoring, Number 1. OECD Principles on Good Laboratory Practice (as revised 1997) ENV/MC/CHEM/(98)17; and European Community (EC) – European

Parliament and council Directive 2004/10/EC (O.J. No. L 50/44,

20/02/2004).

A. Test Material 2,4-Dichlorophenoxyacetic acid

Description: Water Solubility = 569 mg/L

log Kow = 2.81

OECD recommends describing water solubility, melting/boiling point stability in water and light, pKa, Pow or Kow, vapor pressure of test compound, expiration date.

Lot No./Batch No.: Lot No. 2006 2433 8006-USA

Purity: 98.6%

Impurities: None identified

Stability of Compound: The measured 2,4-D concentrations in the replicates across the treatments

ranged from 28.0 to 105% of nominal concentrations. Measured concentrations were highly variable, particularly at lowest test concentrations,

and the study report attributed this to possible microbial metabolism in the

tanks even though there was regular cleaning of the tanks. The CVs for the

replicates in the lowest treatment group ranged from 38-49%. The mean

measured concentrations were 0.245, 3.14, 34.0, and 96.5 mg a.i./L 2,4-D which was 61.3, 78.5, 85.0, and 96.5% of nominal concentrations,

respectively.

Storage Conditions of

Test Chemicals: Not reported

B. Test Organism

Table 2: General Information About the Test Species and Acclimation.

Parameter	Value(s)	Details or Remarks	Guideline Recommendations
Species common name:	Fathead Minnow		EPA recommends fathead minnow (Pimephales promelas).
Species scientific name:	Pimephales promelas		
Species strain (if stated):	not reported		
Were fish obtained from a single laboratory stock?	not reported	Test fish were obtained from the commercial supplier, New England Bioassay, Manchester, CT.	EPA recommends that fish be from a single laboratory stock.
Were acclimation conditions same as definitive test?	Yes	water quality, temperature, and lighting were the same as test conditions	EPA recommends that fish be acclimated under water quality and illumination conditions that are similar to the definitive test.
Acclimation period:	30 days		EPA recommends a minimum two-week acclimation period. Note that the acclimation period is different from the subsequent, in situ pre-exposure phase.

Page 8 of 86

Data Evaluation Record on the Fish Short-Term Reproduction Assay with 2,4-D

Parameter	Value(s)	Details or Remarks	Guideline Recommendations
Details on health:		For 2 weeks out of the 30 days period,	EPA recommends that mortality during the 7
		mortality was reported at a rate of 1.8%	days prior to the pre-exposure phase be less
			than 5% of the culture population. If mortality
			during these 7 days is greater than 10%, EPA
			recommends that the fish be rejected. If
			mortality is between 5-10%, EPA recommends
			that fish be held another 7 days. If mortalities
			greater than 5% occur during this extended
			acclimation period, EPA recommends that the
			fish not be used.
Type of food:	thawed, frozen brine		EPA recommends that fish be fed frozen brine
	shrimp		shrimp twice per day to promote active
			reproduction and maintain body condition.
Source of food:	Brine Shrimp Direct,		
	Ogden, UT		
Frequency of feeding:	at least 2 times/day	1.25-2.0 ml shrimp/replicate vessel	

Page 9 of 86

# Data Evaluation Record on the Fish Short-Term Reproduction Assay with 2,4-D

## EPA MRID Number 48317001

Parameter	Value(s)	Details or Remarks	Guideline Recommendations
Details on feeding:		Feeding start date not reported, but fish	
		were fed during the holding/acclimation,	
		pre-exposure, and exposure periods. Food	
		was withheld 12 hours prior to test	
		termination. Fish food tested free of relevant	
		contaminants.	

Table 3: Fish Selection and Pre-Exposure Performance.

Parameter	Value(s)	Details or Remarks	Guideline Recommendations
Age at test initiation:	5.5 months	6 months at beginning of exposure period	EPA recommends reproductively mature (sexually dimorphic) fish, 4.5 - 6 months old.
Mean weight of males at test initia-tion (if determined):	3.7 g	Mean calculated from all males selected for the test	EPA recommends that a subsample of fish be weighed before the test to estimate the mean weight for each sex.
Range of individual weights (males) at test initiation (if determined):	individual data not pro- vided	range of individual weights were kept within ±20% of the mean	It is recommended that the individual weight of each fish selected for the test be within ±20% of the estimated mean for each sex.
Mean weight of females at test inititation (if determined):	1.9 g	Mean calculated from all females selected for the test	

Page 10 of 86

Data Evaluation Record on the Fish Short-Term Reproduction Assay with 2,4-D

including temperature, photoperiod, feeding, etc., EPA recommends that additional tanks set up at sufficient replicates with the correct sex ratio are EPA recommends that pre-exposure conditions, the beginning of pre-exposure will ensure that EPA recommends a minimum of 14 days. Guideline Recommendations be identical to definitive test conditions. available for the definitive test. The 32 test vessels included 12 test vessels more than the 20 that would be needed for range of individual weights were kept within **Details or Remarks** ±20% of the mean the definitive test. Individual data not Value(s) Not reported Not reported provided 14 days Yes 32 ~ 4 Mean length of females at test initi-Were pre-exposure conditions iden-Mean length of males at test initia-Duration of pre-exposure phase: males) at test initiation (if deter-Range of individual weights (fe-Number of pre-exposure tanks: Number of females per tank: Number of males per tank: tical to the definitive test? Parameter ation (if determined): tion (if determined): mined):

Page 11 of 86

Data Evaluation Record on the Fish Short-Term Reproduction Assay with 2,4-D

Parameter	Value(s)	Details or Remarks	Guideline Recommendations
Pre-exposure fecundity:	>15 eggs/female/reproductive day/replicate	mean = 28.8 eggs/female/day; range from 5.6 to 53.2 eggs/female/day across all replicates	EPA recommends that pre-exposure fecundity in each replicate (tank) selected for use in the definitive test be at least 15 eggs/female/reproductive day/replicate during the 7 days prior to the definitive test.
Number of spawns during pre-ex-posure:	>2 times in 7 days	The top 20 performing spawning groups (>15 eggs/female/reproductive day) were used in the exposure period and were assigned to concentration levels via complete randomized block design.	EPA recommends that spawning occur at least twice in the 7 days prior to the definitive test.
Details on pre-exposure:		Spawning substrate were inspected and the number of eggs laid and eggs found to be infertile were recorded daily	

Page 12 of 86

## Data Evaluation Record on the Fish Short-Term Reproduction Assay with 2,4-D

## EPA MRID Number 48317001

### C. Exposure System

Table 4: Summary of Information on the Exposure System and Test Vessel Characteristics.

Parameter	Value(s)	Details or Remarks	Guideline Recommendations
Type of exposure:	Fflow-through		EPA recommends the use of a flow-through system. As noted in the Corrections and Clarifications document, the use of a static renewal system is not recommended for this assay.
Type of flow-through dilution system:	Continuous flow diluter	A continuous-flow diluter (syringe pump delivery) system delivered stock solution to a mixing chamber where the solution was mixed with dilution water. The test solution was split equally to replicate test vessels.	Intermittent flow proportional diluters or continuous flow serial diluters are recommended.²
Flow-through rate:	45 ± 4.5 mL/min		Recommended flow-through rate is 45 mL/min (2.7 L/hr), or at least 6 total volume exchanges per day.

<sup>&</sup>lt;sup>1</sup> U.S. Environmental Protection Agency (EPA). (2011). Corrections and Clarifications on Technical Aspects of the Test Guidelines for the Endocrine Disruptor Screening Program Tier 1 Assays (OCSPP Test Guideline Series 890). March 3, 2011. Office of Chemical Safety and Pollution Prevention (OCSPP), Washington, D.C. (http://www.epa.gov/endo/pubs/assayvalidation/clarificationdoc.pdf).

Page 13 of 86

<sup>&</sup>lt;sup>2</sup> Additional guidance for aquatic test design is located in OCSPP Guideline 850.1000, Special Considerations for Conducting Aquatic Laboratory Studies.

Data Evaluation Record on the Fish Short-Term Reproduction Assay with 2,4-D

EPA recommends aeration if dissolved oxygen reaches $\leq$ 4.9 mg/L ( $\leq$ 60% saturation).		not reported	Aeration?
	respectively.		
	each treatment level, the %CV for the 0.40 and 100 mg ai./L levels was 0.646% and 0.348%,		
	treatment levels. Of the 4 samples taken at		
	in one vessel at the 0.40 and 100 mg a.i./L		
	Homogeneity of test concentration was analyzed		
10%.			
recommended flow splitting accuracy is within	of each other.		
test system is recommended; 4) The	stock solution-test vessel pairs were within 10%		
completely mixed before introduced into the	splitting accuracy of the concentrations for each		
demonstration that the test solution is	108% to 121%. Within that range, the flow		
not recommended for mixing; 3) A	percent of target concentration ranged from		
recommended but not required; 2) Aeration is	vessels, mixing cells, and stock solutions, the		
through systems: 1) Mixing chamber is	period. For concentrations measured in test		flow-through systems:
Recommended toxicant mixing for flow-	Diluter system was calibrated prior to exposure		Details on toxicant mixing for
Guideline Recommendations	b) Details or Remarks	Value(s)	Parameter

Page 14 of 86

Data Evaluation Record on the Fish Short-Term Reproduction Assay with 2,4-D

Parameter	Value(s)	Details or Remarks	Guideline Recommendations
Source of dilution water:	Natural water	Dilution water originated from Lake Huron which was supplied to the Dow Chemical Company via the City of Midland Water Treatment Plant. The water from the lake was limed and flocculated with ferric chloride and then pumped to the laboratory. Before use in the lab, the water was sand-filtered, pH-adjusted with gaseous CO <sub>2</sub> , carbon-filtered, and UV-irradiated.	EPA recommends natural or reconstituted water: it is recommended that natural water be sterilized with UV and tested for pesticides, heavy metals, and other possible contaminants. OECD accepts any water in which the test species show control survival at least as good as indicated in the test guideline.
Was dilution water analyzed for pesticides, heavy metals, and other contaminants?	Yes	Dilution water is bi-annually analyzed for pesticides, organics, metals and other inorganics.	
Test vessel type/materials:	Glass sealed together with clear silicone adhesive		EPA and OECD recommend that water-contact portions of the system not compromise the study (e.g., all glass vessels or glass vessels with stainless steel frames are acceptable examples).
Test vessel size:	29 cm x 20 cm x 25 cm	approximate fill volume; water depth 13 cm	EPA recommends the use of 18 L test chambers (e.g., 40 x 20 x 20 cm).  EPA recommends 10 L solution per tank.

Page 15 of 86

# Data Evaluation Record on the Fish Short-Term Reproduction Assay with 2,4-D

## EPA MRID Number 48317001

Parameter	Value(s)	Details or Remarks	Guideline Recommendations
Spawning substrate material:	PVC pipe, cut lengthwise		EPA recommends that each tank contain three semi-circular spawning substrates, e.g., aged PVC pipe, 10 - 20 cm in length, split lengthwise.
Spawning substrate size:	9 cm	9 cm long, 10 cm diameter	
Additional details on exposure system:			

Temperature was measured in one test vessel through the exposure period. Temperature, DO, and pH were measured in each test vessel at Day 0 and weekly thereafter during the assay (Table 5). Hardness, alkalinity, and conductivity were measured from a sample in the control and highest exposure level solutions at Day O and weekly thereafter (Table 5).

Table 5: Summary of Water Quality Characteristics in the Test System.

Parameter	Minimum	Maximum	Mean	Measurement Interval	Guideline Recommendations
Temperature (°C)	24.5	25.2	25.0'	Weekly	EPA recommends temperature $25\pm {\it l}^{\circ}C$ ; interreplicate and inter-treatment differentials should not exceed ${\it l}^{\circ}C$ .
	7.02	2.76	7.381	Weekly	EPA recommends pH 6.5 to 9.0.
Dissolved oxygen (mg/L)	6.20	8.29	7.391	Weekly	EPA recommends dissolved oxygen (DO) >4.9 mg/L (>60% air saturation)

Page 16 of 86

Data Evaluation Record on the Fish Short-Term Reproduction Assay with 2,4-D

measured in all test tanks at least weekly and EPA recommends that total organic carbon in EPA recommends that unionized ammonia in continuous temperature monitoring of at least EPA recommends total alkalinity >20 mg/L that hardness and alkalinity be measured in controls and in one tank at the highest test General recommendations for frequency of temperature, pH, and dissolved oxygen be EPA recommends that residual chlorine in concentration at least weekly. In addition, measurements: EPA recommends that Guideline Recommendations the dilution water be  $\leq$ 1 Dg/L. dilution water be <10 Dg/L. dilution water be <2 mg/L. one tank is encouraged. as CaCO3. Measurement Interval Biannually biannually Weekly Weekly weekly Once Mean <200 188.51 36.5 63.51 240 <0.10 mg/L Maximum <2 mg/L (as N) 207 46 89 Minimum 177 62 30 Fotal alkalinity (mg/L as CaCO<sub>3</sub>) Other: Conductivity (µmhos/cm) Hardness [mg/L as CaCO<sub>3</sub>] Total organic carbon (mg/L) Unionized ammonia (Dg/L) Residual chlorine (Dg/L) Other: Ammonia (mg/L) Parameter

1 Means were calculated by the reviewer as the average of the minima and maxima for the ranges provided across control and treated levels

Page 17 of 86

D. Study Design and Additional Experimental Conditions

Table 6: Range-Finding Study Conditions (if Applicable).

Parameter	Value(s)	Details or Remarks	Guideline Recommendations
Was a range-finder conducted?	ON		EPA recommends conducting a range-finder if $96$ -hour LC $_{50}$ data for the fathead minnow are unavailable.
If yes, what was the method for determining the highest test concentration in the rangefinder?	ĄN		EPA recommends that the highest test concentration be selected based on toxicity data for other fish studies or species, if available. Otherwise, either the solubility limit of the test compound or 100 mg/L
Species:	ĄZ		(whichever is lower) is appropriate.
Life stage:	NA		EPA recommends that range-finding tests be performed with fish of similar age and size to those that would be utilized in the test.
Test duration:	NA		EPA recommends a 96-hour exposure.

Page 18 of 86

Data Evaluation Record on the Fish Short-Term Reproduction Assay with 2,4-D

(i)		
iai details:	<b>4</b> 2	EFA recommends conducting a range-finder
		with five test concentrations plus a control
		(six total treatment levels), with four
		females and two males per exposure tank
		(36 fish total). The number of mortalities
		that occur may be used to develop a
		concentration-response curve.
		Based upon the results, the highest
		concentration that does not result in
		increased mortality or signs of overt
		morbidity compared to controls, or 1/3 the
		derived 96-hr LC <sub>50,</sub> may be selected as the
		highest exposure concentration in the 21-
		day test.

Table 7: Definitive Study Conditions.

Parameter	Value(s)	Details or Remarks	Guideline Recommendations
Test duration:	21 days		EPA recommends that the duration of the
			definitive test be 21 days.

Page 19 of 86

Data Evaluation Record on the Fish Short-Term Reproduction Assay with 2,4-D

Guideline Recommendations	EPA recommends that the highest test concentration is either the solubility limit of the test compound, 100 mg/L, or demonstrates adequate evidence of toxicity (e.g., 1/3 the 96-hour LC <sub>50</sub> ), whichever concentration is lowest.		EPA suggests that a concentration separation of between 0.33 (or three-fold) and 0.1 (or tenfold) is scientifically acceptable'.
Details or Remarks	The nominal test concentrations were based on two available 2,4-D studies with the fathead minnow with reported LC <sub>50</sub> values of 320 mg a.i./L (Alexander <i>et al.</i> , 1985) and 133 mg a.i./L (Mayer and Ellersieck, 1986). From these two lethal endpoints, the experimental high concentrations could be estimated at 100 and 40 mg a.i./L, respectively, as specified in the 890.1350 guideline. Because of the uncertainty between the two concentrations, both were selected as high concentrations and 4 concentrations were used for the experimental test.		
Value(s)	Reference studies with fathead minnow	Alexander <i>et al.</i> , 1985; Mayer and Ellersieck, 1986	0, 0.4, 4.0, 40, 100 mg/L
Parameter	Method for selecting the highest test concentration in the definitive test:	Reference study citation (if applicable):	Separation of test concentrations:

Page 20 of 86

Data Evaluation Record on the Fish Short-Term Reproduction Assay with 2,4-D

Parameter	Value(s)	Details or Remarks	Guideline Recommendations
Number of test concentrations:	4		EPA recommends a minimum of 3 concentrations and a control, plus solvent control if appropriate.
Are nominal concentrations adjusted for purity?	Not reported		
Indicate the type of values presented for measured concentrations:	Mean-measured		
Limit of quantification (LOQ):	0.1 mg a.i./L		EPA recommends that for chemical test concentrations below the LOQ, analyses be conducted on the stock solutions.
Level of detection (LOD):	not reported		
Frequency of measurement:	0, 7, 14, and 21 days		It is recommended that test item concentration be measured prior to the addition of fish in all tanks and at least weekly thereafter in two replicates per treatment level.

Page 21 of 86

Data Evaluation Record on the Fish Short-Term Reproduction Assay with 2,4-D

Parameter	Value(s)	Details or Remarks	Guideline Recommendations
Was the randomized complete	Yes		EPA recommends that all fish be randomly
block design used?			assigned to tanks during pre-exposure. Tanks
			fecundity, and the tanks with the highest
			fecundity are randomly assigned to a definitive
			test treatment and block first. Each block
			contains one replicate of each treatment,
			including controls.
Number of replicates in control:	4		EPA recommends 4 replicates.
Number of replicates in solvent	NA		EPA recommends the use of a concurrent
control (if applicable):			solvent control when a solubilizing agent is used.
			EPA recommends 4 replicates.
Number of replicates per test item	4		EPA recommends 4 replicates.
treatment level:			
Number of male fish per replicate at test initiation:	5		EPA recommends 2 males per replicate.
Number of female fish per replicate	4		EPA recommends 4 females per replicate.
at test initiation:			
Was a solvent used?	No		
Was a positive control used?	No		

Page 22 of 86

Data Evaluation Record on the Fish Short-Term Reproduction Assay with 2,4-D

Parameter	Value(s)	Details or Remarks	Guideline Recommendations
Photoperiod:	16 hrs light : 8 hrs dark		EPA recommends photoperiod 16:8 (light:dark).
Light intensity at water's surface:	591 - 770 lux		EPA recommends light intensity 540 – 1080 lux (at water's surface).
Additional details:		Information regarding test solution appearance during the study did not appear to be reported.	

Page 23 of 86

# Data Evaluation Record on the Fish Short-Term Reproduction Assay with 2,4-D

## EPA MRID Number 48317001

Table 8: Summary of Treatment Concentrations in the Fish Short-Term Reproduction Assay with 2,4-D.

	legimoN	Mean Mean			
Treatment ID	Concentration (mg a.i./L)	Concentration (mg a.i./L)	Mean CV (%)	Details or Remarks	Guideline Recommendations
Control (dilution water only)	0.00	<0.10	₹		EPA recommends that test item concentrations be maintained at a coefficient of variation (CV) <>20%.
Treatment 1	0.400	0.245	43.31	61.3% nominal	
Treatment 2	4.00	3.14	13.3	78.5% nominal	
Treatment 3	40.0	34.0	2.0	85.0% nominal	
Treatment 4	100	96.5	4.5	96.5% nominal	

Abbreviations: CV Coefficient of variation.

metabolizing 2,4-D. Though this occurred in all concentrations, the effect was proportionately greater on the lowest nominal concentration. The %CVs were >20% for each time point in this treatment group, though authors report that the variation in concentration was consistent among the replicates and thus the overall exposure 1 Study authors reported a dramatic decrease in measured concentrations between Days 7 and 14 due to a suspected increase in microbial populations capable of was consistent among replicates.

Page 24 of 86

### E. Observations

Biological Endpoints:

Mortality, external abnormalities, abnormal behavior relative to controls, fecundity, fertilization, secondary sex characteristics, full body weight and length, gonadal status (GSI and histology), nuptial tubercle scoring, and vitellogenin concentrations. Observations of appearance were made at test termination, than fish were euthanized, weighed and measured. Blood was collected from the caudal peduncle using heparinized capillary tubes. Fish viscera were removed and fixed in Davidson's solution then preserved in formalin. Gonads were removed, weighed, placed into a plastic tissue cassette and then Davidson's fixative solution. Vitellogenin blood plasma levels were measured with the fathead minnow VTG ELISA kit obtained from Biosense Laboratories, Bergen, Norway. In addition to test samples, each ELISA plate contained 6 calibration standards and 2 non-specific binding assay blanks. Plasma steroid concentrations were not measured.

Were raw (individual) data provided? Yes, with the exception of the re-analysis of female VTG levels and the coloration/appearance of females. Data were also provided in spreadsheet form via email.

EPA recommends that observations of survival, fecundity, fertilization success, secondary sex characteristics, and other clinical signs occur at least daily. At test termination (Day 21), additional observations include body weight and length, nuptial tubercle score, gonadal staging and histopathology, plasma vitellogenin, and plasma sex steroids (testosterone and  $17\beta$ -estradiol, if measured). Gonado-somatic index (GSI) is calculated using a ratio of gonad weight to body weight (gonad weight to the nearest 0.1 mg / body weight in mg x 100) at test termination.

Clinical signs of overt toxicity may include (but are not limited to) hemorrhage, cessation of feeding, and other abnormal behavior.

### II. RESULTS AND DISCUSSION

### A. Results

Only one fish from the 96.5 mg a.i./L treatment group died during the test. No mortalities occurred in the control fish.

Table 9: Adult Fish Survival in Fathead Minnow.

Treatment (mg a.i./L)		Males			Females	
[mean-measured]	n	# Surviving	% Survival	n	# Surviving	% Survival
Negative Control	8	8	100	16	16	100
0.245	8	8	100	16	16	100
3.14	8	8	100	16	16	100
34.0	8	8	100	16	16	100
96.5	8	8	100	16	15	93.8

n = number of individuals per treatment at test initiation.

LOQ=0.10 mg a.i./L  $\,$ 

# Data Evaluation Record on the Fish Short-Term Reproduction Assay with 2,4-D

## EPA MRID Number 48317001

34, and 96.5 mg ai./L). Mean male length ranged from 51.5 mm (negative control) to 54.5 mm(34 mg ai./L), and mean female length ranged from 42.1 mm (3.14 mg Mean male body weight ranged from 3.57 g (negative control) to 3.93 g (34 mg a.i./L), and mean female body weight ranged from 1.59 g (negative control) to 1.62 g (3.14, a.i./L) to 42.5 mm (96.5 mg a.i./L).

Table 10: Size at Test Termination in Fathead Minnow.

			Body	Body Weight					Length	gth		
Treatment (mg a.i./L)		Males			Females			Males			Females	
[mean-measured]	u	Mean (g)	#SD	С	Mean (g)	∓SD	u	Mean (mm)	∓SD	u	Mean (mm)	#SD
Negative Control ( <loq)< td=""><td>4</td><td>3.57</td><td>0.141</td><td>4</td><td>1.59</td><td>0.169</td><td>4</td><td>51.5</td><td>1.54</td><td>4</td><td>42.4</td><td>1.37</td></loq)<>	4	3.57	0.141	4	1.59	0.169	4	51.5	1.54	4	42.4	1.37
0.245	4	3.72	0.308	4	1.64	0.074	4	53.0	1.37	4	42.7	0.56
3.14	4	3.69	0.144	4	1.62	0.173	4	52.2	0.89	4	42.1	0.92
34.0	4	3.93	0.488	4	1.62	0.173	4	54.5	2.20	4	42.4	1.47
96.5	4	3.73	0.408	4	1.62	0.082	4	53.8	2.45	4	42.5	0.78

SD Standard deviation.

n = number of independent replicates per treatment.

LOQ=0.10 mg a.i./L

Page 27 of 86 Version: 22 September 2011

Fecundity ranged from 12.5 to 38.6 eggs/female/reproductive data across all treatment groups, and fertilization success as a percentage of embryos to unfertilized eggs ranged from 91.7 to 99.4%.

Table 11: Fecundity and Fertilization Success in Fathead Minnow.

Treatment	Fecund	lity <sup>1</sup>	Fertilization So	uccess (%) <sup>2</sup>
(mg a.i./L) [mean-measured]	Mean	± SD	Mean	± SD
Negative Control ( <loq))< td=""><td>29.8</td><td>6.4</td><td>96.9</td><td>0.50</td></loq))<>	29.8	6.4	96.9	0.50
0.245	29.0	6.47	96.3	3.08
3.14	25.9	5.27	96.1	1.54
34.0	22.9	8.72	97.7	0.78
96.5	19.7	5.97	95.8	2.68

LOQ=0.10 mg a.i./L

<sup>&</sup>lt;sup>1</sup> Fecundity is calculated as the number of eggs per surviving female per reproductive day per replicate.

<sup>&</sup>lt;sup>2</sup> Fertilization success (%) is calculated as the number of embryos divided by the number of eggs, multiplied by 100.

Male media tubercle scores ranged from 31 (negative control and 0.245 mg a.i./L) to 35 (3.14 mg a.i./L). No tubercles were noted in females.

Table 12: Nuptial Tubercle Score in Fathead Minnow.

Treatment	Ma	lles	Fem	ales <sup>1</sup>
(mg a.i./L) [mean-measured]	n	Median Tubercle Score <sup>2</sup>	n	Median Tubercle Score
Negative Control ( <loq)< td=""><td>4</td><td>31</td><td>4</td><td>0</td></loq)<>	4	31	4	0
0.245	4	31	4	0
3.14	4	35	4	0
34.0	4	33	4	0
96.5	4	34	4	0

n = number of independent replicates per treatment.

### LOQ=0.10 mg a.i./L

<sup>&</sup>lt;sup>1</sup> The study authors reported that no tubercles were observed on female fish.

Mean tubercle scores: 32, 31, 34, 33, and 35 for the negative control and mean-measured 0.245, 3.14, 34.0, and 96.5 mg a.i./L treatment levels, respectively.

Mean male GSI ranged from 1.2 to 1.3%, and mean female GSI ranged from 11.2% (3.15 mg a.i./L) to 13.3% (96.5 mg a.i./L).

Table 13: Gonado-Somatic Index (GSI) in Fathead Minnow.

Treatment		Males			Females	
(mg a.i./L) [mean-measured]	n	Mean GSI <sup>1</sup> (%)	±SD	n	Mean GSI <sup>1</sup> (%)	±SD
Negative Control ( <loq)< td=""><td>4</td><td>1.2</td><td>0.1</td><td>4</td><td>11.9</td><td>1.7</td></loq)<>	4	1.2	0.1	4	11.9	1.7
0.245	4	1.2	0.3	4	12.9	2.3
3.14	4	1.2	0.1	4	11.2	2. 6
34.0	4	1.2	0.2	4	11.8	3.4
96.5	4	1.3	0.1	4	13.3	2.0

n = number of independent replicates per treatment.

LOQ=0.10 mg a.i./L

The study author reports that there were no effects of 2,4-D on the histology or germ cell distribution (staging) of either the testes or ovaries of the test fish. All histopathological findings were considered basal variation unassociated with 2,4-D exposure due to a lack of a concentration-response relationship. Although not concentration-dependent, an increase in the number of female ovaries that were observed as Stage 2 (compared to Stage 3 or 4) was noted by the reviewer for the 2,4-D treatments compared to the negative control. A summary of the reported stage in the ovaries is shown in Table 14b.

The study report did not include observations of decreased proportion of spermatagonia, increased vascular or interstitial proteinaceous fluid, asynchronous gland development, altered proportions of spermatocytes or spermatids, or granulomatous inflammation. There was a single observation of mineralization of the seminiferous tubule/duct (severity grade 1) in one control male. Macrophages (histiocytes; severity grade 1) were observed in female fish in the controls and at all treatment levels.

Gonado-somatic index (%) is calculated as gonad weight (to the nearest 0.1 mg) / body weight (mg) x 100.

Table 14a: Gonadal Staging in Fathead Minnow.

Treatment	Ма	les	Fem	ales
(mg a.i./L) [mean-measured]	n	Median Stage <sup>1</sup>	n	Median Stage <sup>2</sup>
Negative Control ( <loq)< td=""><td>8</td><td>3</td><td>4</td><td>4</td></loq)<>	8	3	4	4
0.245	8	4	4	4
3.14	8	3	4	3
34.0	4	3	4	3
96.5	4	3	4	4

n = total number of animals per treatment upon which observations were made at test termination.

### LOQ=0.10 mg a.i./L

Table 14b. Ovary staging (from study report)

Ovary	Control	0.245 mg a.i./L	3.14 mg a.i./L	34 mg a.i./L	96.5 mg a.i./L
Examined	16	16	16	16	15
Stage 2	0	3	4	6	3
Stage 3	6	2	6	3	4
Stage 4	10	11	6	7	8

The guideline recommends the following gonadal staging scale for male fathead minnow: O=undeveloped, 1=early spermatogenic, 2=mid-spermatogenic, 3=late spermatogenic, 4=spent.

The guideline recommends the following gonadal staging scale for female fathead minnow: O=undeveloped, 1=early development, 2=mid-development, 3=late development, 4=late development/hydrated, 5=post-ovulatory.

Table 15: Gonadal Histopathology in Male Fathead Minnow.

					Diagnost	ic Obser	rvations <sup>1</sup>				
Treatment (mg a.i./L) [mean- measured]	Severity	Pro	creased portion of rmatogonia		esence of estis-Ova	semir	ization of niferous bule	Т	ncreased esticular generation	Нур	rstitial Cell pertrophy/ perplasia
		n	Incidence	n	Incidence	n	Incidence	n	Incidence	n	Incidence
Negative	0	8	8	8	8	8	7	8	8	8	8
Control	1	8	0	8	0	8	1	8	0	8	0
( <loq)< td=""><td>2</td><td>8</td><td>0</td><td>8</td><td>0</td><td>8</td><td>0</td><td>8</td><td>0</td><td>8</td><td>0</td></loq)<>	2	8	0	8	0	8	0	8	0	8	0
	3	8	0	8	0	8	0	8	0	8	0
	4	8	0	8	0	8	0	8	0	8	0
0.245	0	8	8	8	8	8	8	8	8	8	8
	1	8	0	8	0	8	0	8	0	8	0
	2	8	0	8	0	8	0	8	0	8	0
	3	8	0	8	0	8	0	8	0	8	0
	4	8	0	8	0	8	0	8	0	8	0
3.14	0	8	7	8	8	8	8	8	8	8	8
	1	8	1	8	0	8	0	8	0	8	0
	2	8	0	8	0	8	0	8	0	8	0
	3	8	0	8	0	8	0	8	0	8	0
	4	8	0	8 0		8	0	8	0	8	0
34.0	0	8	7	8	8	8	8	8	8	8	8
	1	8	0	8	0	8	0	8	0	8	0
	2	8	0	8	0	8	0	8	0	8	0
	3	8	1	8	0	8	0	8	0	8	0
	4	8	0	8	0	8	0	8	0	8	0
96.5	0	8	8	8	8	8	8	8	8	8	8
	1	8	0	8	0	8	0	8	0	8	0
	2	8	0	8	0	8	0	8	0	8	0
	3	8	0	8	0	8	0	8	0	8	0
	4	8	0	8	0	8	0	8	0	8	0

LOQ=0.10 mg a.i./L

Gonadal histopathology diagnostic observations are graded 0 - 4 based on severity: 0=Not remarkable, 1=Minimal, 2=Mild, 3=Moderate, 4=Severe. See Appendix E of the test guideline for reference.

Data Evaluation Record on the Fish Short-Term Reproduction Assay with 2,4-D

Table 16: Additional Gonadal Histopathology Observations in Male Fathead Minnow.

					Additional Diagnostic Observations <sup>2</sup>	nostic (	Observations <sup>2</sup>				
Treatment		De	Decreased	Increa	Increased Vascular	Asyı	Asynchronous	Altere	Altered Proportions	Gran	Granulomatoris
(mg a.i./L)	Severity	Pro	Proportion of	or	or Interstitial		Gonad	of S	of Spermatocytes	i ii	Inflammation
[mean-measured]		Sper	Spermatogonia	Proteil	Proteinaceous Fluid	Dev	Development	or	or Spermatids		
		<sup>L</sup> u	Incidence	n <sup>1</sup>	Incidence	r_ L	Incidence	Ľ	Incidence	Ē	Incidence
Negative Control	0	NA	AN	NA	NA	A A	ΑN	N A	ΑN	A A	AN
( <foø)< td=""><td>1</td><td>NA</td><td>NA</td><td>NA</td><td>NA</td><td>NA</td><td>NA</td><td>NA</td><td>NA</td><td>NA</td><td>NA</td></foø)<>	1	NA	NA	NA	NA	NA	NA	NA	NA	NA	NA
	2	NA	NA	NA	NA	NA	NA	NA	VΝ	NA	NA
	3	NA	NA	NA	NA	NA	NA	NA	NA	NA	NA
	4	NA	NA	NA	NA	NA	NA	NA	NA	NA	NA
0.245	0	NA	NA	NA	NA	NA	NA	NA	NA	NA	NA
	1	NA	NA	NA	NA	NA	NA	NA	NA	NA	NA
	2	NA	NA	NA	NA	NA	NA	NA	NA	NA	NA
	3	NA	NA	NA	NA	NA	NA	NA	NA	NA	NA
	4	AN	NA	NA	NA	NA	NA	NA	VΑ	NA	NA

Page 34 of 86

Version: 22 September 2011

Data Evaluation Record on the Fish Short-Term Reproduction Assay with 2,4-D

EPA MRID Number 48317001

			ons Granulomatous es Inflammation			Granul Inflam n <sup>1</sup>	Granul Inflam NA NA	Granul Inflam NA NA	Granul Inflam NA NA NA NA	Granul Inflam NA NA NA NA	Granul Inflam NA NA NA NA	Granul Inflam NA NA NA NA NA NA	Granul Inflam NA NA NA NA NA NA NA NA NA NA	Granul Inflam NA NA N
	Altered Proportions	of Spermatocytes	or Spermatids	n <sup>1</sup> Incidence	NA	NA	NA	NA	NA	NA	NA	NA	NA	NA
Observations <sup>2</sup>	Asynchronous	Gonad	Development	Incidence	NA	AN	NA	NA	NA	NA	NA	NA	NA	AN
nostic	As		۵	n	NA	A	NA	NA	A	NA	NA	NA	ΑΝ	A
Additional Diagnostic Observations <sup>2</sup>	Increased Vascular	or Interstitial	Proteinaceous Fluid	Incidence	ΥN	ΥN	ΝΑ	ΝΑ	ΥN	ΝΑ	ΝΑ	ΝΑ	ΥN	ΥN
	Increa	or	Proteil	<sup>_</sup> u	NA	AN	NA	NA	AN	NA	NA	NA	AN	AN
	Decreased	Proportion of	Spermatogonia	Incidence	NA	AN	NA	NA	AN	NA	NA	NA	NA	ΑN
	De	Pro	Sper	r L	NA	AN	NA	NA	AN	NA	NA	NA	AN	AN
		Severity	Severity		0	1	2	3	4	0	1	2	3	4
	Treatment	(mg a.i./L)	[mean-measured]		3.14					34.0				

Page 35 of 86

Version: 22 September 2011

Data Evaluation Record on the Fish Short-Term Reproduction Assay with 2,4-D

					Additional Diagnostic Observations <sup>2</sup>	nostic (	)bservations <sup>2</sup>				
Treatment		De	Decreased	Increa	Increased Vascular	Asyr	Asynchronous	Altere	Altered Proportions	30.	
(mg a.i./L)		Prog	Proportion of	o	or Interstitial		Gonad	of Sp	of Spermatocytes	5 5	indiidiomiation
[mean-measured]	Severity	Speri	Spermatogonia	Protein	Proteinaceous Fluid	Dev	Development	o	or Spermatids		
		n¹	Incidence	n¹	Incidence	n	Incidence	n <sup>1</sup>	Incidence	ιu	Incidence
96.5	0	AN	ΑN	NA	AN	A N	AN	₹ Z	ΑN	NA	NA
	1	NA	ΑN	NA	ΑN	AN	NA	AN	ΑN	VΝ	NA
	7	NA	NA	NA	NA	NA	NA	NA	NA	VΝ	NA
	3	NA	NA	NA	NA	NA	NA	NA	NA	NA	NA
	4	NA	NA	NA	NA	AN	NA	AN	NA	NA	NA

Abbreviations: Not applicable.

LOQ=0.10 mg a.i./L

<sup>1</sup> Number of individual fish observed

Gonadal histopathology diagnostic observations are graded 0 - 4 based on severity: 0=Not remarkable, 1=Minimal, 2=Mild, 3=Moderate, 4=Severe. See

Appendix E of the test guideline for reference.

Page 36 of 86

Version: 22 September 2011

Table 17: Gonadal Histopathology in Female Fathead Minnow.

			Diagno	ostic Obs	servations <sup>1</sup>		
Treatment (mg a.i./L) [mean-measured]	Severity		ased Oocyte Atresia	H	follicular Cell yperplasia/ ypertrophy		creased Yolk Formation
		n	Incidence	n	Incidence	n	Incidence
Negative Control	0	16	9	16	16	16	16
( <loq)< td=""><td>1</td><td>16</td><td>7</td><td>16</td><td>0</td><td>16</td><td>0</td></loq)<>	1	16	7	16	0	16	0
	2	16	0	16	0	16	0
	3	16	0	16	0	16	0
	4	16	0	16	0	16	0
0.245	0	16	8	16	16	16	16
	1	16	5	16	0	16	0
	2	16	3	16	0	16	0
	3	16	0	16	0	16	0
	4	16	0	16	0	16	0
3.14	0	16	10	16.	16	16	16
	1	16	3	16	0	16	0
	2	16	2	16	0	16	0
	3	16	1	16	0	16	0
	4	16	0	16	0	16	0
34.0	0	16	7	16	16	16	16
	1	16	5	16	0	16	0
	2	16	3	16	0	16	0
	3	16	1	16	0	16	0
	4	16	0	16	0	16	0
96.5	0	15	11	15	15	15	15
	1	15	3	15	0	15	0
	2	15	1	15	0	15	0
	3	15	0	15	0	15	0
	4	15	0	15	0	15	0

LOQ=0.10 mg a.i./L

Table 18: Additional Gonadal Histopathology Observations in Female Fathead Minnow.

Treatment			А	ddition	al Diagnostic (	Observa	ations <sup>1</sup>		
(mg a.i./L) [mean-	Severity	Interst	itial Fibrosis		g Debris in Oviduct		anulomatous		eased Post- tory Follicles
measured]		n	Incidence	n	Incidence	n	Incidence	n	Incidence
Negative	0	16	16	16	16	16	16	16	16
Control	1	16	0	16	0	16	0	16	0
( <loq)< td=""><td>2</td><td>16</td><td>0</td><td>16</td><td>0</td><td>16</td><td>0</td><td>16</td><td>0</td></loq)<>	2	16	0	16	0	16	0	16	0
	3	16	0	16	0	16	0	16	0
	4	16	0	16	0	16	0	16	0
0.245	0	16	16	16	16	16	16	16	16
	1	16	0	16	0	16	0	16	0
	2	16	0	16	0	16	0	16	0
	3	16	0	16	0	16	0	16	0
	4	16	0	16	0	16	0	16	0
3.14	0	16	16	16	16	16	15	16	16
	1	16	0	16	0	16	1	16	0
	2	16	0	16	0	16	0	16	0
	3	16	0	16	0	16	0	16	0
	4	16	0	16	0	16	0	16	0
34.0	0	16	16	16	16	16	16	16	16
	1	16	0	16	0	16	0	16	0
	2	16	0	16	0	16	0	16	0
	3	16	0	16	0	16	0	16	0
	4	16	0	16	0	16	0	16	0

Gonadal histopathology diagnostic observations are graded 0 - 4 based on severity: 0=Not remarkable, 1=Minimal, 2=Mild, 3=Moderate, 4=Severe. See Appendix E of the test guideline for reference.

Treatment			А	ddition	al Diagnostic (	Observ	ations <sup>1</sup>		
(mg a.i./L) [mean-	Severity	Interst	itial Fibrosis		Debris in		ranulomatous nflammation		eased Post- tory Follicles
measured]		n	Incidence	n	Incidence	n	Incidence	n	Incidence
96.5	0	15	15	15	15	15	15	15	15
	1	15	0	15	0	15	0	15	0
	2	15	0	15	0	15	0	15	0
	3	15	0	15	0	15	0	15	0
	4	15	0	15	0	15	0	15	0

Gonadal histopathology diagnostic observations are graded 0 - 4 based on severity: 0=Not remarkable, 1=Minimal, 2=Mild, 3=Moderate, 4=Severe. See Appendix E of the test guideline for reference.

LOQ=0.10 mg a.i./L

Mean male plasma vitellogenin ranged from 0.909 ng/mL (34 mg a.i./L) to 1.77 ng/mL (0.245 mg a.i./L). Mean female plasma vitellogenin ranged from  $19.8 \times 10^6$  ng/mL (3.14 mg a.i./L) to  $55.1 \times 10^6$  ng/mL (0.245 mg a.i./L).

Table 19: Plasma Vitellogenin in Fathead Minnow.

			Plasma Vite	ellogenin (	(VTG)	
Treatment (mg a.i./L)		Males			Females	
[mean-measured]	n	Mean (ng/mL plasma)	±SD	n	Mean (ng/mL plasma)	±SD
Negative Control ( <loq)< td=""><td>4</td><td>1.34</td><td>1.73</td><td>4</td><td>27.6×10<sup>6</sup></td><td>21.5×10<sup>6</sup></td></loq)<>	4	1.34	1.73	4	27.6×10 <sup>6</sup>	21.5×10 <sup>6</sup>
0.245	4	1.77	1.25	4	55.1×10 <sup>6</sup>	47.5×10 <sup>6</sup>
3.14	4	0.930	0.777	4	19.8×10 <sup>6</sup>	7.69×10 <sup>6</sup>
34.0	4	0.909	0.591	4	21.0×10 <sup>6</sup>	14.4×10 <sup>6</sup>
96.5	4	1.29	0.983	4	31.6×10 <sup>6</sup>	12.7×10 <sup>6</sup>

SD Standard deviation.

n = number of independent replicates per treatment.

 $LOQ=0.10 \ mg \ a.i./L$ 

## Data Evaluation Record on the Fish Short-Term Reproduction Assay with 2,4-D

### EPA MRID Number 48317001

Plasma testosterone and plasma  $17\beta$ -estradiol in male and females were not measured (Table 20).

Table 20: Plasma Sex Steroids in Fathead Minnow (Pimephales promelas). Not measured.

		Д.	lasma Te	Plasma Testosterone (T)	E			Plas	Plasma 17β-estradiol (E2)	estradiol	(E2)	
Treatment		Males			Females			Males			Females	
(mg a.i./L) [mean-measured]	د	Mean (ng/mL plasma)	#SD	c	Mean (ng/mL plasma)	∓SD	۵	Mean (ng/mL plasma)	#SD	c	Mean (ng/mL plasma)	#SD
Negative Control ( <loq)< td=""><td>₹ Z</td><td>Ω</td><td>Q N</td><td>∢ Z</td><td>QN</td><td>ΩN</td><td>₹ Z</td><td>Q.</td><td>Q.</td><td>₹ Z</td><td>Q</td><td>ON O</td></loq)<>	₹ Z	Ω	Q N	∢ Z	QN	ΩN	₹ Z	Q.	Q.	₹ Z	Q	ON O
0.245	Ž	QN	Q.	Ϋ́	QN	ND	ΑN	Q.	ND	ΑΝ	QN	ND
3.14	Ž	QN	Q.	Υ Υ	QN	ND	ΑN	9	ND	Ϋ́	Q	ND
34.0	A A	QN	ND	۸N	QΝ	QN	AN	QN	ND	ΑN	QN	ND
96.5	A A	QN	ND	AN	QN	ND	AN	ND	ND	NA	QN	ND

Abbreviations: NA Not applicable. ND Not determined. SD Standard deviation.

LOQ=0.10 mg a.i./L

quiescence, and feeding abstinence. It was reported that one female fish in the 0.245 mg a.i./L treatment group lost an eye during the exposure No abnormal behavior was observed among control or treatment fish, such as hyperventilation, loss of equilibrium, uncoordinated swimming, atypical

Page 41 of 86 Version: 22 September 2011

## Data Evaluation Record on the Fish Short-Term Reproduction Assay with 2,4-D

### EPA MRID Number 48317001

period which could have been from breeding activity and/or colliding with breeding substrate. Another female fish in the 3.14 mg a.i./L treatment treatment group and one fish in the control group. Vertical banding was observed in both control and treatment group female test fish, but no group exhibited ascites (i.e., accumulation of fluid in the coelomic cavity). Scoliosis, or bent tail, was displayed by one fish in the 96.5 mg a.i./L treatment-related effects were observed. The raw data for this banding was not reported.

Table 21: Secondary Sex Characteristics and Clinical Signs in Fathead Minnow at Test Termination.

Treatment		Seconda	ıry Sex Chara	Secondary Sex Characteristics and Clinical Signs		
(mg a.i./L)	Males			Females		
[mean-measured]	Туре	n	Incidence	Туре	u	Incidence
Negative Control ( <loq)< td=""><td>None</td><td>ΝΑ</td><td>ΥN</td><td>None</td><td>NA</td><td>NA</td></loq)<>	None	ΝΑ	ΥN	None	NA	NA
0.245	None	NA	NA	None	NA	NA
3.14	None	ΥN	ΨN	None	₹ Z	Ϋ́
34.0	None	NA	ΝΑ	None	ΑN	ΑN
96.5	None	NA	Ϋ́Ν	None	AN	NA

Abbreviations: NA Not applicable.

LOQ=0.10 mg a.i./L

Page 42 of 86

B. Study Author's Analysis and Conclusions

The study author reported the following statistical methodology:

The appropriate units of statistical analyses were the measures of central tendency from the replicate test vessels. The statistical significance of all tests was judged at the 0.05 significance level, with the exception of the Shapiro-Wilk test, which was judged at the 0.01 significance level. All biological response data, apart

from mortality, were analyzed and reported separately by sex.

Endpoints were statistically evaluated with the Jonckheere-Terpstra test (Hollander and Wolfe, 1973) in a

step down manner if the measures of central tendency for that endpoint were consistent with a monotone

concentration-response. If the endpoint response was not consistent with a monotone concentration-

response, the data were assessed for normality using the Shapiro-Wilk test and variance homogeneity using

Levene's test. Where non-normality or variance heterogeneity was observed, normalizing and/or variance

stabilizing transformations were applied. If the data were non monotonic and normally distributed with

homogeneous variances, then a significant treatment effect was determined using the one way ANOVA

followed by Dunnett's test. Where non-normality or variance heterogeneity was observed normalizing and/or

variance stabilizing transformations was applied. If the data were normally distributed with homogeneous

variances then a significant treatment effect was determined using the one way ANOVA followed by Dunnett's

test. If the data were normally distributed with heterogeneous variance, the Mann-Whitney-Wilcoxon U test

was used. If no normalizing transformation was found, the Mann-Whitney-Wilcoxon U test using a

Bonferroni-Holm adjustment was used.

Significant mortality was assessed, if necessary, using the Cochran-Armitage Linear Trend Test where the

data was consistent with a monotonic concentration response, and otherwise from Fisher's Exact test with

a Bonferroni-Holm adjustment. A treatment effect for tubercle score was determined using the multiquantal

Jonckheere-Terpstra test.

Female plasma VTG was reanalyzed by the study authors in an effort to demonstrate proficiency in the

method as outlined in the 890.1350 guideline. Results of the reanalysis are presented in Appendix II to

this DER.

The enzyme-linked immunosorbent assay (ELISA) test was used to analyze the vitellogenin (VTG) levels

in male and female test fish blood plasma. An ELISA reanalysis of VTG concentrations in female fish

demonstrated the ability of using the ELISA in accordance with the methods recommended by the 890.1350

guideline; insufficient male plasma was available for reanalysis. For the primary test, the correlation

coefficient (R2) for the calibration curves were all >0.99 (meeting 890.1350 standards); the concentrations

for each standard fell between 82.1 and 115% of nominal concentrations (meeting 890.1350 standards);

and the average absorbance measured for the non-specific binding assay blanks was less than 0.06

absorbance units (highest recommended level by 890.1350 guidelines) although one plate produced an

average absorbance unit value of 0.067.

For the reanalysis of female VTG samples, the R2 were all >0.99 and the concentrations for each standard

fell between 70-120% of nominal (850.1350 recommended range) with the exception of one standard with

a recovery of 144%. Two aliquots of each sample dilution were analyzed on the same ELISA plate (well

duplicates) to verify the co-efficient of variation. It was reported that 74 of the 79 samples did not differ by

more than 20%, and all were less than 31%, and thus these results indicated acceptable method precision.

The actual response in the well-duplicate samples was greater than the primary concentration analysis. The

study author indicated that it has been shown that concentrations of VTG may increase when re-analyzed

following a freeze-thaw cycle (Brodeur, et al., 2006).

C. Reviewer's Analysis and Conclusions

Statistical Methods: Statistical analyses were recalculated by the reviewer following decision logic provided

in test guidelines. Continuous data were tested for normality using Shapiro-Wilks test and for homogeneity

of variance using Levene's test. Data that demonstrated a monotonic trend were analyzed using Jonckheere-

Terpstra according the 890.1350 guideline flowchart and using Williams test according to EFED policy.

Data that met the assumptions of normality and homogeneity of variance were then analyzed using Dunnett's

test if the data did not demonstrate a monotonic concentration response. Data that failed the assumption of

normality or homogeneity of variance were analyzed using the Mann-Whitney-Wilcoxon test. Mortality was

assessed using the Fisher's exact test with a Bonferroni-Holm adjustment.

Conclusions: Fecundity was significantly decreased in the highest concentration tested, 96.5 mg a.i./L

compared to the control. There was no significant difference for fertilization success (# embryos/# eggs

x 100), male tubercle score, male and female condition index (whole body length/weight), male and female GSI, or male and female VTG levels. No signs of abnormal behavior, appearance, or secondary sex characteristics in the treatment groups relative to the control. Additionally, the study author reports that there were no effects of 2,4-D on the histology or germ cell distribution (staging) of either the testes or ovaries of the test fish. While not concentration-dependent, an increase in the number of female ovaries that were observed as Stage 2 (compared to Stage 3 or 4) was reported for the 2,4-D treatments compared to the control; the number of Stage 2 ovaries was 0, 3, 4, 6, and 3 for the control, 0.245, 3.14, 34.0 and 96.5 mg a.i./L treatment group, respectively.

# Data Evaluation Record on the Fish Short-Term Reproduction Assay with 2,4-D

## **EPA MRID Number 48317001**

Table 22: Reproductive and HPG Endpoints<sup>1,2,3</sup> for Male Fathead Minnow in the FSTRA with 2,4-D.

					Gonadal						
Treatment	Tubercle Score	Score	GSI	<u></u>	Staging and	Plasn	Plasma VTG	Plasma T	аТ	Plasma E2	L E2
(mg a.i./L)					Histo.						
[mean-measured]	20:10 M	2	#:0 %	2	Effect?	#!'0 %	٥	#!U 70	2	#!()	2
	Mediali	ı.		D.	(Yes/No)	, C	2.	% <u>C</u>	<b>1</b>	% <b>DIII.</b>	Ъ
Negative Control	21	< 2	\ <u>\</u>	2	Š	< 2	\ \ \ \ \ \ \ \ \ \ \ \ \ \ \ \ \ \ \	Š	2	\ <u>\</u>	\ \ \ \ \ \ \ \ \ \ \ \ \ \ \ \ \ \ \
( <fog)< td=""><td><u>-</u></td><td>Ž</td><td>Í Ž</td><td><u> </u></td><td><u> </u></td><td><u> </u></td><td><u> </u></td><td>Ĭ Ž</td><td><u> </u></td><td><u>(</u> Z</td><td>ζ Z</td></fog)<>	<u>-</u>	Ž	Í Ž	<u> </u>	<u> </u>	<u> </u>	<u> </u>	Ĭ Ž	<u> </u>	<u>(</u> Z	ζ Z
0.245	31	0.995	2.29	866.0	ON	32.1	0.954	VΝ	ΨN	ΑN	NA
3.14	35	0.925	-1.25	1.000	ON	-30.3	0.962	ΝΑ	ΨN	ΑN	NA
34.0	33	0.978	-1.25	1.000	ON	-32.1	0.954	ΝΑ	ΨN	ΑN	NA
96.5	34	0.748	8.32	0.817	No	-3.6	1.000	NA	NA	NA	NA
Statistical Test	Dunnett's <sup>4</sup>	tt's⁴	Dunnett's	ett's	AN	Dun	Dunnett's	AN		AN	

NA Not applicable

LOQ=0.10 mg a.i./L

1 Unless otherwise indicated, effects and percent (%) differences are reported based on comparison to the negative (clean water) control. Conclusions regarding histopathology may be heavily weighted by the expert opinion of a board-certified pathologist.

<sup>2</sup> Unless otherwise specified, effects are considered statistically significant at p<0.05.

<sup>3</sup> For percent (%) difference, positive values indicate an increase relative to the negative control, and negative values indicate a decrease relative to the negative control.

<sup>4</sup> Mean tubercle scores: 32, 31, 34, 33, and 35 for the negative control and mean-measured 0.245, 3.14, 34.0, and 96.5 mg a.i./L treatment levels, respectively.

Page 46 of 86

# Data Evaluation Record on the Fish Short-Term Reproduction Assay with 2,4-D

## **EPA MRID Number 48317001**

Table 23: Reproductive and HPG Endpoints<sup>1,2,3</sup> for Female Fathead Minnow in the FSTRA with 2,4-D.

Treatment	Fecu	Fecundity	Fert. Si	Fert. Success	Tubercle Score	Score	ő	ISSI	Gonadal Staging and	Plasms	Plasma VTG	Plasma T	na T	Plasma E2	E2
(mg a.i./L)									Histo.						
[mean-	%	1	%	1	9	1	%	1	Effect?	%	1	%	1	ä 2	•
liedsuled]	Diff.	ď	Diff.	Д	Median	Д	Diff.	ď	(Yes/No)	Diff.	Д	Diff.	д	% DIII.	d
Negative Control	2	4	4	4	c	4	<u> </u>	4	4	2	4	4	4	\ - 4	<u> </u>
( <loq)< td=""><td>Į Z</td><td>Į Ž</td><td>Į Ž</td><td>¥ Z</td><td><b>&gt;</b></td><td>¥ Ž</td><td>¥ Z</td><td>ď Z</td><td>۲ 2</td><td>₹ Z</td><td>¥ Ž</td><td><u> </u></td><td><u> </u></td><td><u> </u></td><td><u> </u></td></loq)<>	Į Z	Į Ž	Į Ž	¥ Z	<b>&gt;</b>	¥ Ž	¥ Z	ď Z	۲ 2	₹ Z	¥ Ž	<u> </u>	<u> </u>	<u> </u>	<u> </u>
0.245	-2.6	0.332	-0.64	0.975	0	AN	8.21	0.945	No	9.66	0.494	NA	NA	ĄN	Ą
3.14	-13	0.210	-0.83	0.943	0	ΑN	-6.11	0.980	No	-28.4	0.889	NA	NA	ΑN	Ą
34.0	-23	0.088	22.0	0.954	0	NA	-1.13	1.00	No	-23.5	8/9.0	NA	NA	AN	AN
96.5	-34	0.031	-1.16	0.839	0	NA	11.5	0.847	No	14.6	0.494	NA	NA	AN	AN
	Jonck	Jonckheere-		,44	4			,	<u> </u>	1		2	<	2	
Statistical Test	Ter	Terpstra	Dunneu s	ell s	¥ Z		IUNO	Dunnett s	<u> </u>	Mann-	Mann-vynitney	ZY.	1	YN.	

Abbreviations:  $^{\text{Diff.}}$  Difference.  $^{\text{E2}}$  17 $\beta$ -estradiol.  $^{\text{Fert}}$  Fertilization.  $^{\text{GSI}}$  Gonado-Somatic Index.  $^{\text{Histo}}$  Histopathology.

NA Not applicable. T Testosterone. VTG Vitellogenin.

LOQ=0.10 mg a.i./L

1 Unless otherwise indicated, effects and percent (%) differences are reported based on comparison to the negative (clean water) control. Conclusions regarding histopathology may be heavily weighted by the expert opinion of a board-certified pathologist.

Page 47 of 86

 $<sup>^2\,</sup>$  Unless otherwise specified, effects are considered statistically significant at p<0.05.

<sup>&</sup>lt;sup>3</sup> For percent (%) difference, positive values indicate an increase relative to the negative control, and negative values indicate a decrease relative to the negative control.

Table 24: Growth Endpoints<sup>1,2,3</sup> in the Fish Short-Term Reproduction Assay (FSTRA) with 2,4-D.

Treatment		Body \	Weight			Le	ength	
(mg a.i./L)	Ma	les	Fem	ales	Mal	es	Fer	nales
[mean-measured]	% Diff.	р	% Diff.	р	% Diff.	р	% Diff.	р
Negative Control ( <loq)< td=""><td>NA</td><td>NA</td><td>NA</td><td>NA</td><td>NA</td><td>NA</td><td>NA</td><td>NA</td></loq)<>	NA	NA	NA	NA	NA	NA	NA	NA
0.245	4.20	0.914	2.67	0.979	2.96	0.579	0.59	0.992
3.14	3.36	0.958	2.04	0.992	1.26	0.959	-0.71	0.984
34.0	10.01	0.380	2.04	0.992	5.82	0.097	-0.12	1.000
96.5	4.48	0.895	1.57	0.997	4.32	0.273	0.18	1.000
Statistical Test	Dunn	ett's	Duni	nett's	Dunn	ett's	Dur	nett's

LOQ=0.10 mg a.i./L

### E. Study Deficiencies

Several deviations are listed in Section I. Materials and Methods of this DER. One test performance criterion was not met; the CV for the mean-measured concentration of the lowest treatment group was 43% which is greater than the guideline criterion of <20%. These deviations did not impact the interpretation of this study.

### F. Reviewer's Comments

Two replicate spawning groups with performance ranking within the top 20 breeding groups were inadvertently excluded from assignment to the exposure period test vessels. The study author reported that this deviation did not impact the experimental test because the 20 spawning groups used for the exposure period met the minimum criteria as outlined in the 890.1350 guideline.

<sup>&</sup>lt;sup>1</sup> Unless otherwise indicated, percent (%) differences are reported based on comparison to the negative (clean water) control.

<sup>&</sup>lt;sup>2</sup> Unless otherwise specified, effects are considered statistically significant at p<0.05.

For percent (%) difference, positive values indicate an increase relative to the negative control, and negative values indicate a decrease relative to the negative control.

Ammonia and chlorine levels were measured in dilution water biannually. Ammonia levels were detected at

0.24 mg a.i./L, exceeding guideline recommendations of  $\leq 1$  µg/L. Additionally, OCSPP guidelines

recommend chlorine levels <10  $\mu g/L$  but chlorine in the dilution water was only quantified at a level <200

μg/L (chlorine was not detectable at 200 μg/L). Organic carbon levels in dilution water were not measured.

The study author's rationale for the decrease in test concentrations in the test vessels but not the stock

solutions was the marked increase in growth and waste production from Day 7 to Day 14 which increased

microbial populations. Test material microbial biodegradation had a greater percentage effect in the two

lower treatment groups than the two highest treatment groups. The study author reported that biodegradation

is known to occur under the conditions in the current study (Sinton et al., 1986). It was noted that the

variability in 2,4-D concentrations in the stock solutions might have occurred because the large stock

volumes were not being continuously mixed during their use.

The reviewer reports the gonad histopathological results as they were reported by the study author (Table

6). The stage and severity grading of the gonads was conducted in accordance with the procedures

described in the 890.1350 guidelines. Gonad sections were examined under a light microscope by an

American College of Veterinary Pathologists board-certified veterinary pathologist.

Page 49 of 86

### III. REFERENCES

- Alexander, H.C., Gersich, F.M., and Mayes, M.A. (1985). Acute Toxicity of Four Phenoxy Herbicides to Aquatic Organisms. *Bull. Env. Contam. Toxicol.* 35, 314-321.
- Brodeur, J.C., Woodburn K.B., Zhang F., Bartels M.J., and Klecka G.M. (2006). Plasma sampling and freezing procedures influence vitellogenin measurements by enzyme-linked immunoassay in the fathead minnow \*Pimephales promelas). Env. Toxicol. Chem. 25, 337-348.
- Hansch, C., Leo, A., and Hoekman, D. (1995). Exploring QSAR Hydrophobic, Electronic and Steric Constants, ACS, Washington, D.C.
- Hopkins, D.L. (1987). 2,4-Dichlorophenoxyacetic Acid: Determination of the water solubility. Report of Toxicology & Environmental Research and Consulting, The Dow Chemical Company, Midland Michigan.
- Lezotte F., Van Hoven R., and Nixon, W. (2002). Determination of Water Solubility of 2,4-D Acid by the Shake Flask Method: Lab Project Number: 467C-103: 467/121801/7840/SUB467. Unpublished study prepared by Wildlife International, Ltd. MRID No. 45692501.
- Mayer, F. M. and Ellersieck, M.R. (1986). Manual of Acute Toxicity: Interpretation and Data Base for 410 Chemicals and 66 Species of Freshwater Animals. Resource Publication 160. U.S. Department of the Interior, Fish and Wildlife Service. Washington D.C.
- Mayes, M. A., Gorzinski, S. J., Potter, R. B., and Richardson, C. H. (1990). 2,4-Dichlorophenoxyacetic acid: Evaluation of the Toxicity to Early Life Stages of the Fathead Minnow, Pimephales promelas. Report of Toxicology & Environmental Research and Consulting, The Dow Chemical Company, Midland, Michigan.
- Sinton, G. L., Fan, L. T., Erickson, L. E., and Lee, S. M. (1986). Biodegradation of 2,4-D and related xenobiotic compounds. Enzyme Microb. Technol. 8, 395-403.
- Wo, C. (2008). Certificate of Analysis for 2,4-D Acid, Lot# 2006 2433 8006-USA. Study number 26625, PSL Reference No. 081210-2H. Eurofms|Product Safety Laboratories, Dayton, New Jersey.

Data Evaluation Record on the Fish Short-Term Reproduction Assay with 2,4-D

APPENDIX I. OUTPUT OF REVIEWER'S STATISTICAL VERIFICATION:

Endpoint	Monotonic?	Parametric?	890.1350	EFED	Comments
Survival			Fisher's Exact Test	Fisher's Exact Test	
Overall:	No	ΑΝ	n.s p>0.05	n.s p>0.05	No significant effect on overall survival
Female body weight	o <sub>N</sub>	Yes	Dunnett's:	Dunnett's:	890.1350 and EFED same conclusions, no
			n.s. p>0.05	n.s. p>0.05	effect
Male body weight	oN	Yes	Dunnett's:	Dunnett's:	890.1350 and EFED same conclusions, no
			n.s. p>0.05	n.s. p>0.05	effect
Female body length	oN	Yes	Dunnett's:	Dunnett's:	890.1350 and EFED same conclusions, no
			n.s. p>0.05	n.s. p>0.05	effect
Male body length	No	Yes	Dunnett's:	Dunnett's:	890.1350 and EFED same conclusions, no
			n.s. p>0.05	n.s. p>0.05	effect
Female VTG	Yes	Yes	Mann Whitney;	Mann Whitney;	890.1350 and EFED same conclusions, no
			n.s. p>0.05	n.s. p>0.05	effect
Male VTG	Yes	No	Dunnett's:	Dunnett's:	890.1350 and EFED same conclusions, no
			n.s. p>0.05	n.s. p>0.05	effect
Female GSI	Yes	Yes	Dunnett's:	Dunnett's:	890.1350 and EFED same conclusions, no
			n.s. p>0.05	n.s. p>0.05	effect
Male GSI	No	Yes	Dunnett's:	Dunnett's:	890.1350 and EFED same conclusions, no
			n.s. p>0.05	n.s. p>0.05	effect
Female tubercle	NA	NA	NA	NA	No score
score					

Page 51 of 86

Data Evaluation Record on the Fish Short-Term Reproduction Assay with 2,4-D

EPA MRID Number 48317001

Endpoint	Monotonic?	Parametric?	890.1350	EFED	Comments
Male tubercle score	Yes	Yes	Dunnett's:	Dunnett's:	890.1350 and EFED same conclusions, no
			n.s. p>0.05	n.s. p>0.05	effect
Fecundity	No	Yes	Jonckheere-Terpstra:	Williams;	890.1350 and EFED same conclusions, effect
			dose 3 p=0.011	dose 3 p=0.031	at dose 3
Fertility	Yes	No	Dunnett's:	Dunnett's:	890.1350 and EFED same conclusions, no
			n.s. p>0.05	n.s. p>0.05	effect
F testosterone	NA	NA	NA	NA	NA
M testosterone	NA	NA	NA	NA	NA
F estradiol	NA	NA	NA	NA	NA
M estradiol	NA	NA	NA	NA	NA

Page 52 of 86

### Appendix II: Reanalysis of Female Plasma VTG.

Table: Plasma Vitellogenin in Fathead Minnow (Results of Reanalysis of Female Samples).

		Plasma Vitellogenin (VTG)	
Treatment (mg a.i./L)		Females	
[measured]	n	Mean (ng/mL plasma)	±SD
Negative Control	4	59.8x10 <sup>6</sup>	26.6x10 <sup>6</sup>
0.245	4	74.8x10 <sup>6</sup>	31.0x10 <sup>6</sup>
3.14	4	55.7x10 <sup>6</sup>	18.8x10 <sup>6</sup>
34.0	4	58.7x10 <sup>6</sup>	33.7x10 <sup>6</sup>
96.5	4	92.5x10 <sup>6</sup>	27.6x10 <sup>6</sup>

SD Standard deviation.

n = number of independent replicates per treatment

### Appendix III: Reviewer's Statistical Output.

	_			
Fi	cher	' a	Exact	Teat

		NUMBER OF				
IDENTIFICATION	ALIVE	DEAD	TOTAL ANIMALS			
CONTROL	16	0	16			
0.245	16	0	16			
TOTAL	32	0	32			

\_\_\_\_\_\_

Critical Fisher's value (16,16,16) (alpha=0.05) is 11.0. b value is 16. Since b is greater than 11.0 there is no significant difference between CONTROL and TREATMENT at the 0.05 level.

### Fisher's Exact Test

		NUMBER OF			
IDENTIFICATION	ALIVE	DEAD	TOTAL ANIMALS		
CONTROL	16	0	16		
3.14	16	0	16		
TOTAL	32	0	32		

\_\_\_\_\_\_

Critical Fisher's value (16,16,16) (alpha=0.05) is 11.0. b value is 16. Since b is greater than 11.0 there is no significant difference between CONTROL and TREATMENT at the 0.05 level.

Fisher's Exact Test

Page 54 of 86

		NUMBER OF				
IDENTIFICATION	ALIVE	DEAD	TOTAL ANIMALS			
CONTROL	16	0	16			
34.0	16	0	16			
TOTAL	32	0	32			

Critical Fisher's value (16,16,16) (alpha=0.05) is 11.0. b value is 16. Since b is greater than 11.0 there is no significant difference between CONTROL and TREATMENT at the 0.05 level.

Fisher's Exact Test

		NU:	HBER OF
IDENTIFICATION	ALIVE	DEAD	TOTAL ANIMALS
CONTROL	16	0	16
96.5	15	1	16
TOTAL	31	1	32

Critical Fisher's value (16,16,16) (alpha=0.05) is 11.0. b value is 15. Since b is greater than 11.0 there is no significant difference between CONTROL and TREATMENT at the 0.05 level.

Summary of Fisher's Exact Tests

GROUP	IDENTIFICATION	NUMBER EXPOSED	NUMBER DEAD	SIG 0.05
	CONTROL	16	0	
1	0.245	16	0	
2	3.14	16	0	
3	34.0	16	0	
4	96.5	16	1	

test for fish screen study - TEST DATA 2 4 D ANALYSIS RESULTS FOR VARIABLE VAR01 ( F body weight (g) )

TESTS OF ASSUMPTIONS FOR PARAMETRIC ANALYSIS Shapiro-Wilks test for Normality of Residuals -- alpha-level=0.01 Levenes test for homogeneity of variance(absolute residuals) -- alpha-level=0.05 Use parametric analyses if neither test rejected, otherwise non-parametric analyses.

Shapiro-Wilks Shapiro-Wilks Levenes Levenes Conclusion

	Stat 952		value 0.390	Test Sta 2.884		E PARAMETRIC	TESTS
*****	****	*****	*****	*****	*****	*****	*****
BASIC SU	MMARY	STATIST	CICS				
Level	N	Mean	StdDev	StdErr	Coef of Var	95% Conf.I	nterval
Ctrl	4	1.59	0.17	0.08	10.51	1.33,	1.86
Dose1	4	1.64	0.07	0.04	4.48	1.52,	1.75
Dose2	4	1.63	0.17	0.09	10.65	1.35,	1.90
Dose3	4	1.63	0.17	0.09	10.52	1.35,	1.90
Dose4	4	1.62	0.08	0.04	5.11	1.49,	1.75
Level		Median	Min	Max	%of Control(mean	s) %Reduct	ion(means)
Ctrl		1.53	1.47	1.84			
Dose1		1.61	1.58	1.74	102.67	-2.6	7
Dose2		1.64	1.43	1.79	102.04	-2.0	4
Dose3		1.63	1.45	1.80	102.04	-2.0	4
Dose4		1.62	1.52	1.72	101.57	-1.5	7

Analysis of Variance (ANOVA) - overall F-test

Numerator df Denominator df F-stat P-value 4 15 0.05 0.994

Dunnett - testing each trt mean signif. different than control Williams - test assumes dose-response relationship, testing negative trend Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC

Level Mean Dunnett Isotonic Williams Tukey p-values
p-value mean p-value Dose1 Dose2 Dose3 Dose4 Dose5

Page 56 of 86

						EPA	MRID Numb	Del 46317C	01
Ctrl	1.59		1.62						
Dose1	1.64	0.979	1.62	0.692					
Dose2	1.63	0.992	1.62	0.727	1.000				
Dose3	1.63	0.992	1.62	0.746	1.000	1.000			
Dose4	1.62	0.997	1.62	0.750	1.000	1.000	1.000	•	
****	*****	*****	*****	*****	*****	*****	*****	*****	****
_	_	ANALYSES		lpha-level			sts		
				among tre		roups			
Deg	-	Freedom	TestStat						
	4		0.38	0.98	4				
lannWhit	t - test	ing each t	ert median	signif. d	ifferent	from co	ntrol		
onckhee	ere - te	st assumes	s dose-res	ponse rela	tionship	, testin	g negativ	ve trend	
Level	Medi	an	MannWhi	t p-value		Jonckh	eere p-va	alue	
Ctrl	1.	53							
Dose1	1.	61		0.343			0.877		
Dose2	1.			1.000			0.722		
Dose3	1.	63		0.889			0.593		
Dose4	1.	62		0.676			0.579		
DECREAS	SING TRE	ND TEST SU	JMMARY	LOWEST C	ONCENTRA:	rion sig	NIF. LESS	S THAN C	ONTRO
Willi	iams						(no sign		
Jonck	cheere				>highe	st dose	(no sign	. differ	ences
*****	*****	*****	*****	*****	*****	*****	*****	*****	***
	RIC ANAL		-	-level=0.0 overall F		l tests			
	_			f F-sta		P-value			
Null	4	1!		0.05		0.994			
	1	Δ,	,	0.03	`	3.221			
				gnif. diff					
			_	nse relati	_	_			
ukey -	two-sid	ed tests,	all possil	ble compar	isons, no	ot used	for NOEC	or LOEC	
evel	Mean	Dunnett	Isotonic	Williams			Tukey p-v	values	
		p-value	mean	p-value	Dose1	Dose2	Dose3	Dose4	Dos
Ctrl	-1.59		-1.59						
Dose1	-1.64	0.979	-1.63	0.441	-		-		•
Dose2	-1.63	0.992	-1.63	0.470	1.000	•	•	•	•
Dose3	-1.63	0.992	-1.63	0.487	1.000	1.000	•	•	•
Dose3	-1.62	0.997	-1.63	0.497	1.000	1.000	1.000	•	•
20201	1.02	0.557	1.03	0.157	1.000	1.000	1.000	•	•
*****	*****	* * * * * * * * *	*****	*****	****	*****	*****	*****	***
				lpha-level					

MannWhit - testing each trt median signif. different from control

0.38

Kruskal-Wallis test - equality among treatment groups

Degrees of Freedom TestStat P-value

0.984

	Jonckheere -	test	assumes	dose-response	relationship,	testing	INCREASING	trend
--	--------------	------	---------	---------------	---------------	---------	------------	-------

Level	Median	MannWhit p-value	Jonckheere p-value
Ctrl	-1.53		
Dose1	-1.61	0.343	0.123
Dose2	-1.64	1.000	0.278
Dose3	-1.63	0.889	0.407
Dose4	-1.62	0.676	0.421

INCREASING TREND TEST SUMMARY CONTROL

LOWEST CONCENTRATION SIGNIF. GREATER THAN

Williams Jonckheere >highest dose (no sign. differences)
>highest dose (no sign. differences)

test for fish screen study - TEST DATA 2 4 D ANALYSIS RESULTS FOR VARIABLE VAR02 ( M body weight (g) )

TESTS OF ASSUMPTIONS FOR PARAMETRIC ANALYSIS

Shapiro-Wilks test for Normality of Residuals -- alpha-level=0.01

Levenes test for homogeneity of variance(absolute residuals) -- alpha-level=0.05 Use parametric analyses if neither test rejected, otherwise non-parametric analyses.

Shapiro-Wilks	Shapiro-Wilks	Levenes	Levenes	Conclusion
Test Stat	P-value	Test Stat	P-value	
0.970	0.762	2.364	0.100	USE PARAMETRIC TESTS

\*\*\*\*\*\*\*\*\*\*\*\*\*\*\*\*\*\*\*

BASIC ST	JMMAR	Y STATIST	TICS				
Level	N	Mean	StdDev	StdErr	Coef of Var	95% Conf.Ir	nterval
Ctrl	4	3.57	0.14	0.07	3.99	3.34,	3.80
Dose1	4	3.72	0.31	0.15	8.32	3.23,	4.21
Dose2	4	3.69	0.14	0.07	3.91	3.46,	3.92
Dose3	4	3.93	0.49	0.24	12.44	3.15,	4.70
Dose4	4	3.73	0.41	0.20	10.91	3.08,	4.38
Level		Median	Min	Max	<pre>%of Control(means)</pre>	) %Reducti	ion(means)
Ctrl		3.53	3.45	3.77			
Dose1		3.69	3.38	4.12	104.20	-4.20	)
Dose2		3.65	3.57	3.90	103.36	-3.36	5
Dose3		3.85	3.45	4.57	110.01	-10.01	L
Dose4		3.75	3.25	4.17	104.48	-4.48	3

\*

PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
Analysis of Variance (ANOVA) - overall F-test
Numerator df Denominator df F-stat P-value
4 15 0.61 0.660

Dunnett - testing each trt mean signif. different than control Williams - test assumes dose-response relationship, testing negative trend Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC

Level Mean Dunnett Isotonic Williams Tukey p-values

					EPA MRID Number 48317001					
		p-value	mean	p-value	Dose1	Dose2	Dose3	Dose4	Dose5	
Ctrl	3.57		3.73							
Dose1	3.72	0.914	3.73	0.826						
Dose2	3.69	0.958	3.73	0.855	1.000					
Dose3	3.93	0.380	3.73	0.870	0.895	0.842				
Dose4	3.73	0.895	3.73	0.879	1.000	1.000	0.911			
				*****				*****	* * * *	
Krus	kal-Wal	ANALYSES lis test - Freedom		lpha-level: among trea P-value 0.76	atment g e		ests			
				signif. da ponse rela				ve trend		
Level Ctrl	Median 3.53		MannWhi		Jonckheere p-value					
Dose1	3.	69		0.678			0.718			
Dose2	3.	65		0.283			0.830			
Dose3	3.	85		0.413			0.920			
Dose4	3.	75		0.678			0.841			
Willi	_	ND TEST SU	JMMARY	LOWEST CO	>highe	st dose	NIF. LES (no sign (no sign	. differ	ences)	
******	*****	*****	*****	*****	*****	*****	*****	*****	***	
Anal		Variance df Deno	(ANOVA) - ominator d		-test t :	P-value				
	4	15	)	0.61		0.660				
Williams	- test	assumes d	lose-respo	gnif. diffense relationse relationse relationse compari	onship,	testing	INCREASI			
Level	Mean	Dunnett p-value	Isotonic mean		Dose1	Dose2	Tukey p-	values Dose4	Dose5	

Level	Mean	Dunnett	Isotonic	Williams	Tukey p-values					
		p-value	mean	p-value	Dose1	Dose2	Dose3	Dose4	Dose5	
Ctrl	-3.57		-3.57	•						
Dose1	-3.72	0.914	-3.71	0.340						
Dose2	-3.69	0.958	-3.71	0.364	1.000	•	•	•		
Dose3	-3.93	0.380	-3.83	0.186	0.895	0.842				
Dose4	-3.73	0.895	-3.83	0.190	1.000	1.000	0.911			

\*\*\*\*\*\*\*\*\*\*\*\*\*\*\*\*\*\*\*

NON-PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests Kruskal-Wallis test - equality among treatment groups Degrees of Freedom TestStat P-value 4 1.83 0.766

Page 59 of 86

MannWhit - testing each trt median signif. different from control Jonckheere - test assumes dose-response relationship, testing INCREASING trend

Level	Median	MannWhit p-value	Jonckheere p-value
Ctrl	-3.53		•
Dose1	-3.69	0.678	0.282
Dose2	-3.65	0.283	0.170
Dose3	-3.85	0.413	0.080
Dose4	-3.75	0.678	0.159

CONTROL

INCREASING TREND TEST SUMMARY LOWEST CONCENTRATION SIGNIF. GREATER THAN

Williams

>highest dose (no sign. differences) Jonckheere >highest dose (no sign. differences)

test for fish screen study - TEST DATA 2 4 D ANALYSIS RESULTS FOR VARIABLE VAR03 ( F body length (mm) )

TESTS OF ASSUMPTIONS FOR PARAMETRIC ANALYSIS Shapiro-Wilks test for Normality of Residuals -- alpha-level=0.01 Levenes test for homogeneity of variance(absolute residuals) -- alpha-level=0.05 Use parametric analyses if neither test rejected, otherwise non-parametric analyses.

Shapiro-Wilks	Shapiro-Wilks	Levenes	Levenes	Conclusion
Test Stat	P-value	Test Stat	P-value	
0.939	0.225	1.586	0.229	USE PARAMETRIC TESTS

\*

BASIC ST	JMMARY	STATIS'	TICS				
Level	N	Mean	StdDev	StdErr	Coef of Var	95% Conf.I	nterval
Ctrl	4	42.43	1.33	0.67	3.14	40.30,	44.55
Dose1	4	42.68	0.57	0.28	1.33	41.77,	43.58
Dose2	4	42.13	0.92	0.46	2.19	40.66,	43.59
Dose3	4	42.38	1.49	0.74	3.51	40.01,	44.74
Dose4	4	42.50	0.81	0.41	1.91	41.21,	43.79
Level		Median	Min	Max	%of Control(means)	%Reduct	ion(means)
Ctrl		42.20	41.10	44.20		•	
Dose1		42.65	42.10	43.30	100.59	-0.5	9
Dose2		42.45	40.80	42.80	99.29	0.7	1
Dose3		42.25	40.90	44.10	99.88	0.1	.2
Dose4		42.85	41.30	43.00	100.18	-0.1	.8

\*

PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests Analysis of Variance (ANOVA) - overall F-test Numerator df Denominator df F-stat P-value 0.14 15

Dunnett - testing each trt mean signif. different than control Williams - test assumes dose-response relationship, testing negative trend Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC

#### EPA MRID Number 48317001

						EPA	MRID Numi	per 4831/0	01
Level	Mean		Isotonic				Tukey p-		
		p-value	mean	p-value	Dose1	Dose2	Dose3	Dose4	Dose5
Ctrl	42.43		42.55						
Dose1	42.68	0.992	42.55	0.651	•		•		
Dose2			42.33	0.565	0.948				
Dose3	42.38	1.000	42.33	0.583	0.994	0.997			
Dose4	42.50	1.000	42.33	0.595	0.999	0.987	1.000		
*****	*****	*****	*****	*****	*****	******	*****	*****	***
NON-PAR	AMETRIC	ANALYSES	- use a	alpha-level	=0.05 fo	r all te	ests		
				among tre					
		Freedom				-			
	4		1.03	0.90	5				
26 771 1		1				_			
				n signif. d sponse rela				ve trend	
OOHOHHIC	CIC CC	ese assanic	.b dobe iei	sponse rera	CIOHDHIP	, ccbcii	ig negaci	ve erena	
Level	Medi		MannWhi	it p-value		Jonckl	neere p-va	alue	
Ctrl	42.								
Dose1				0.678			0.718		
Dose2				1.000			0.442		
Dose3				0.889			0.463		
Dose4	42.	. 85		0.676			0.566		
DECREA	SING TRE	END TEST S	HIMMARY	LOWEST C	ONCENTRA	TION SIG	NIF. LES	S THAN C	ONTROL
Will		IND IEST C	011111111	LOWED1 C			(no sign		
	kheere						(no sign		
				******			******	*****	***
				a-level=0.0		l tests			
	_			overall F		_			
Nu	merator			df F-sta		P-value			
	4	1	.5	0.14		0.966			
Dunnett	- testi	ing each t	rt mean si	gnif. diff	erent tha	an contr	rol		
				onse relati				NG trend	
				ble compar					
Level	Mean	Dunnett	Isotonic				Tukey p-		
		p-value	mean	p-value	Dose1	Dose2	Dose3	Dose4	Dose5
Ctrl	-42.43		-42.40						
Dosel		0.992	-42.40	0.597	•	•	•	•	•
	-42.13		-42.40	0.632	0.948	•	•	•	•
	-42.38		-42.40	0.651	0.994	0.997	•	•	•
	-42.50		-42.50	0.605	0.999		1.000		
				******				*****	***
		ANALYSES		alpha-level			ests		
				among tre		roups			
De	_	Freedom	TestStat						
	4		1.03	0.90	5				

Page 61 of 86

MannWhit - testing each trt median signif. different from control Jonckheere - test assumes dose-response relationship, testing INCREASING trend

Level	Median	MannWhit p-value	Jonckheere p-value
Ctrl	-42.20		•
Dose1	-42.65	0.678	0.282
Dose2	-42.45	1.000	0.558
Dose3	-42.25	0.889	0.537
Dose4	-42.85	0.676	0.434

INCREASING TREND TEST SUMMARY CONTROL

LOWEST CONCENTRATION SIGNIF. GREATER THAN

Williams

>highest dose (no sign. differences) >highest dose (no sign. differences) Jonckheere

test for fish screen study - TEST DATA 2 4 D ANALYSIS RESULTS FOR VARIABLE VAR04 ( M body length (mm) )

TESTS OF ASSUMPTIONS FOR PARAMETRIC ANALYSIS

Shapiro-Wilks test for Normality of Residuals -- alpha-level=0.01 Levenes test for homogeneity of variance(absolute residuals) -- alpha-level=0.05 Use parametric analyses if neither test rejected, otherwise non-parametric analyses.

Shapiro-Wilks	Shapiro-Wilks	Levenes	Levenes	Conclusion
Test Stat	P-value	Test Stat	P-value	
0.938	0.222	1.295	0.316	USE PARAMETRIC TESTS

^ ^ ^ ^ ^ ^ ^ ^ ^ ^ /		^ ^ ^ ^ ^ ^ ^ ^ ^ ^ ^ ^			^ ^ ^ ^ ^ ^ ^ ^ ^ ^ ^ ^ ^ ^ ^ ^ ^ ^ ^ ^			
BASIC SU	JMMAR	Y STATIST	TICS					
Level	N	Mean	StdDev	StdErr	Coef of Var	95% Conf.	Interval	
Ctrl	4	51.55	1.54	0.77	2.98	49.10,	54.00	
Dose1	4	53.08	1.35	0.67	2.53	50.93,	55.22	
Dose2	4	52.20	0.91	0.45	1.74	50.75,	53.65	
Dose3	4	54.55	2.21	1.11	4.06	51.03,	58.07	
Dose4	4	53.78	2.46	1.23	4.58	49.85,	57.70	
Level		Median	Min	Max	%of Control(means)	) %Reduct	cion(means)	
Ctrl		51.55	49.70	53.40				
Dose1		52.60	52.10	55.00	102.96	-2.9	96	
Dose2		52.00	51.40	53.40	101.26	-1.2	26	
Dose3		54.85	52.00	56.50	105.82	-5.8	32	
Dose4		54.50	50.30	55.80	104.32	-4.3	32	

\*

PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests Analysis of Variance (ANOVA) - overall F-test Numerator df Denominator df F-stat P-value 4 15 1.80 0.181

Dunnett - testing each trt mean signif. different than control Williams - test assumes dose-response relationship, testing negative trend Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC

Level	Mean	Dunnett p-value	Isotonic mean	Williams p-value	Dose1	Dose2	Tukey p-v Dose3	values Dose4	Dose5
Ctrl Dose1 Dose2 Dose3	51.55 53.08 52.20 54.55			0.927 0.944 0.953			· ·		
NON-PARA	AMETRIC	ANALYSES	- use a	0.958 ************************************	****** =0.05 fo	****** r all te	* * * * * * * *	• * * * * * * * * *	***
		Tis test Freedom		r among trea P-value 0.19	e	roups			
Jonckhe	ere - te	est assume	s dose-res	signif. d		, testing	g negativ		
Level	Medi	.an	MannWhi	t p-value		Jonckh	eere p-va	alue	

Median	MannWhit p-value	Jonckheere p-value
51.55	•	
52.60	0.235	0.926
52.00	0.580	0.745
54.85	0.125	0.966
54.50	0.235	0.975
	51.55 52.60 52.00 54.85	51.55       .         52.60       0.235         52.00       0.580         54.85       0.125

DECREASING T	REND	TEST	SUMMARY	LOWEST	CONCENTRATIO	N SIC	SNIF.	LESS	THAN	CONTROL
Williams					>highest	dose	(no	sign.	diffe	erences)
Jonckheere	2				>highest	dose	(no	sign.	diffe	erences)

Analysis of Variance (ANOVA) - overall F-test
Numerator df Denominator df F-stat P-value
4 15 1.80 0.181

Dunnett - testing each trt mean signif. different than control Williams - test assumes dose-response relationship, testing INCREASING trend Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC

Level	Mean	Dunnett	Isotonic	Williams		-	Tukey p-v	values	
		p-value	mean	p-value	Dose1	Dose2	Dose3	Dose4	Dose5
Ctrl	-51.55		-51.55	•					•
Dose1	-53.08	0.579	-52.64	0.241					
Dose2	-52.20	0.959	-52.64	0.257	0.955				
Dose3	-54.55	0.097	-54.16	0.035	0.769	0.378	•	•	
Dose4	-53.78	0.273	-54.16	0.036	0.980	0.726	0.971	•	

\*\*\*\*\*\*\*\*\*\*\*\*\*\*\*\*\*\*\*

NON-PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests Kruskal-Wallis test - equality among treatment groups Degrees of Freedom TestStat P-value

Page 63 of 86

4 6.12 0.190 MannWhit - testing each trt median signif. different from control Jonckheere - test assumes dose-response relationship, testing INCREASING trend Median Level MannWhit p-value Jonckheere p-value -51.55 Ctrl 0.074 0.235 Dose1 -52.60 Dose2 -52.00 0.580 0.255 Dose3 -54.85 0.125 0.034 -54.50 0.235 0.025 Dose4 INCREASING TREND TEST SUMMARY LOWEST CONCENTRATION SIGNIF. GREATER THAN CONTROL Williams Dose3 Jonckheere Dose3 test for fish screen study - TEST DATA 2 4 D ANALYSIS RESULTS FOR VARIABLE VAR05 (F vitellogenin (ng/mL)) TESTS OF ASSUMPTIONS FOR PARAMETRIC ANALYSIS Shapiro-Wilks test for Normality of Residuals -- alpha-level=0.01 Levenes test for homogeneity of variance(absolute residuals) -- alpha-level=0.05 Use parametric analyses if neither test rejected, otherwise non-parametric analyses. Shapiro-Wilks Shapiro-Wilks Levenes Levenes Conclusion Test Stat P-value Test Stat P-value 0.172 0.019 0.932 4.120 USE NON-PARAMETRIC TESTS \* BASIC SUMMARY STATISTICS Level N Mean StdDev StdErr Coef of Var 95% Conf.Interval Ctrl 4 27595000 21470025 10735012 77.80 -6568601,61758601 -2.044E7,1.3061E8 86.17 Dose1 4 55085000 47464483 23732241 Dose2 4 19765000 7687585 3843793 38.89 7532336,31997664 Dose3 4 21102500 14406782 7203391 68.27 -1821905,44026905 Dose4 4 31632500 12736128 6368064 40.26 11366478,51898522 Max %of Control(means) Level Median Min %Reduction(means) 20850000 10090000 58590000 Ctrl Dosel 47025000 9630000 1.1666E8 199.62 -99.62 71.63 Dose2 20845000 9690000 27680000 28.37 Dose3 14660000 12450000 42640000 76.47 23.53 Dose4 29195000 18870000 49270000 114.63 -14.63\* PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests Analysis of Variance (ANOVA) - overall F-test Numerator df Denominator df F-stat P-value 4 15 1.30 0.315 Dunnett - testing each trt mean signif. different than control

Williams - test assumes dose-response relationship, testing negative trend

Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOE	Tukey -	two-sided	tests,	all	possible	comparisons,	not	used	for	NOEC	or	LOEC
---	---------	-----------	--------	-----	----------	--------------	-----	------	-----	------	----	------

Level Mean		Dunnett	Isotonic	Williams	Tukey p-values					
		p-value	mean	p-value	Dose1	Dose2	Dose3	Dose4	Dose5	
Ctrl	27595000		41340000							
Dose1	L55085000	0.373	41340000	0.852						
Dose2	219765000	0.976	24166667	0.532	0.316	•		•		
Dose3	321102500	0.988	24166667	0.550	0.350	1.000				
Dose4	131632500	0.998	24166667	0.561	0.682	0.960	0.974			

\*

NON-PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests

Kruskal-Wallis test - equality among treatment groups

Degrees of Freedom TestStat P-value

3.43 0.489

MannWhit - testing each trt median signif. different from control Jonckheere - test assumes dose-response relationship, testing negative trend

Median	MannWhit p-value	Jonckheere p-value
20850000		
47025000	0.494	0.807
20845000	0.889	0.442
14660000	0.678	0.257
29195000	0.494	0.655
		20850000

DECREASING TREND TEST SUMMARY

Williams

Jonckheere

LOWEST CONCENTRATION SIGNIF. LESS THAN CONTROL

>highest dose (no sign. differences)

>highest dose (no sign. differences)

\*

PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
Analysis of Variance (ANOVA) - overall F-test
Numerator df Denominator df F-stat P-value
4 15 1.30 0.315

Dunnett - testing each trt mean signif. different than control Williams - test assumes dose-response relationship, testing INCREASING trend Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC

Level	Mean	Dunnett	Isotonic	Williams	Tukey p-values					
		p-value	mean	p-value	Dose1	Dose2	Dose3	Dose4	Dose5	
Ctrl	-2.76E7	•	-2.76E7	•				•		
Dose1	-5.509E7	0.373	-3.19E7	0.479						
Dose2	-1.977E7	0.976	-3.19E7	0.510	0.316					
Dose3	-2.11E7	0.988	-3.19E7	0.527	0.350	1.000	•			
Dose4	-3.163E7	0.998	-3.19E7	0.538	0.682	0.960	0.974	•		

\*

NON-PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests Kruskal-Wallis test - equality among treatment groups

Page 65 of 86

Degrees of Freedom TestStat P-value 3.43 0.489 MannWhit - testing each trt median signif. different from control Jonckheere - test assumes dose-response relationship, testing INCREASING trend Level MannWhit p-value Jonckheere p-value Median Ctrl -2.085E7 Dose1 -4.703E70.494 0.193 Dose2 -2.085E7 0.889 0.558 Dose3 -1.466E7 0.678 0.743 Dose4 -2.92E7 0.494 0.345 INCREASING TREND TEST SUMMARY LOWEST CONCENTRATION SIGNIF. GREATER THAN CONTROL Williams >highest dose (no sign. differences) Jonckheere >highest dose (no sign. differences) test for fish screen study - TEST DATA 2 4 D ANALYSIS RESULTS FOR VARIABLE VAR06 ( M vitellogenin (ng/mL) ) TESTS OF ASSUMPTIONS FOR PARAMETRIC ANALYSIS Shapiro-Wilks test for Normality of Residuals -- alpha-level=0.01 Levenes test for homogeneity of variance(absolute residuals) -- alpha-level=0.05 Use parametric analyses if neither test rejected, otherwise non-parametric analyses. Shapiro-Wilks Shapiro-Wilks Levenes Levenes Conclusion Test Stat P-value Test Stat P-value 0.048 0.135 USE PARAMETRIC TESTS 0.904 2.073 \* BASIC SUMMARY STATISTICS Level N Mean StdDev StdErr Coef of Var 95% Conf.Interval Ctrl 4 1337.00 1735.75 867.87 129.82 -1424.96, 4098.96 Dosel 4 1766.25 1249.84 624.92 70.76 -222.52, 3755.02 Dose2 4 932.00 779.35 389.67 83.62 -308.11, 2172.11 Dose3 4 907.75 589.11 294.55 64.90 -29.65, 1845.15 Dose4 4 1288.50 981.96 490.98 76.21 -274.01, 2851.01 95% Conf.Interval Level Median Min Max %of Control(means) %Reduction(means) 

 Median
 Min
 Max
 % of Control(r

 492.00
 424.00
 3940.00
 .

 1769.00
 447.00
 3080.00
 132.11

 807.50
 153.00
 1960.00
 69.71

 Ctrl Dose1 -32.11Dose2 30.29 380.00 1750.00 67.89 Dose3 750.50 32.11 1177.00 2440.00 Dose4 360.00 96.37 \* PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests Analysis of Variance (ANOVA) - overall F-test Numerator df Denominator df F-stat P-value 15 0.38 0.819

Dunnett - testing each trt mean signif. different than control

Williams - test assumes dose-response relationship, testing negative trend Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC

Level	Mean	Dunnett	Isotonic	Williams	Tukey p-values					
		p-value	mean	p-value	Dose1	Dose2	Dose3	Dose4	Dose5	
Ctrl	1337.00		1551.63							
Dose1	1766.25	0.954	1551.63	0.691						
Dose2	932.00	0.962	1042.75	0.456	0.835	•				
Dose3	907.75	0.954	1042.75	0.472	0.821	1.000				
Dose4	1288.50	1.000	1042.75	0.482	0.974	0.991	0.989			

\*

NON-PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests

Kruskal-Wallis test - equality among treatment groups

Degrees of Freedom TestStat P-value

1.44 0.837

MannWhit - testing each trt median signif. different from control Jonckheere - test assumes dose-response relationship, testing negative trend

Level	Median	MannWhit p-value	Jonckheere p-value
Ctrl	492.00		
Dose1	1769.00	0.678	0.718
Dose2	807.50	0.889	0.500
Dose3	750.50	0.889	0.390
Dose4	1177.00	0.889	0.447

DECREASING TREND TEST SUMMARY

Williams

Jonckheere

LOWEST CONCENTRATION SIGNIF. LESS THAN CONTROL

\*highest dose (no sign. differences)

>highest dose (no sign. differences)

\*

PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
Analysis of Variance (ANOVA) - overall F-test
Numerator df Denominator df F-stat P-value
4 15 0.38 0.819

Dunnett - testing each trt mean signif. different than control Williams - test assumes dose-response relationship, testing INCREASING trend Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC

Level Mean		Dunnett	Isotonic	Williams	Tukey p-values					
		p-value	mean	p-value	Dose1	Dose2	Dose3	Dose4	Dose5	
Ctrl -	1337.00	•	-1235.75	•		•	•	•	•	
Dosel-	1766.25	0.954	-1235.75	0.636		•				
Dose2	-932.00	0.962	-1235.75	0.671	0.835	•	•			
Dose3	-907.75	0.954	-1235.75	0.690	0.821	1.000				
Dose4-	1288.50	1.000	-1288.50	0.674	0.974	0.991	0.989			

\*

NON-PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests

```
Kruskal-Wallis test - equality among treatment groups
 Degrees of Freedom TestStat P-value
                       1.44
                               0.837
```

MannWhit - testing each trt median signif. different from control Jonckheere - test assumes dose-response relationship, testing INCREASING trend

Level	Median	MannWhit p-value	Jonckheere p-value
Ctrl	-492.00		
Dose1	-1769.00	0.678	0.282
Dose2	-807.50	0.889	0.500
Dose3	-750.50	0.889	0.610
Dose4	-1177.00	0.889	0.553

INCREASING TREND TEST SUMMARY LOWEST CONCENTRATION SIGNIF. GREATER THAN

Williams Jonckheere >highest dose (no sign. differences) >highest dose (no sign. differences)

test for fish screen study - TEST DATA 2 4 D ANALYSIS RESULTS FOR VARIABLE VAR07 ( F GSI )

TESTS OF ASSUMPTIONS FOR PARAMETRIC ANALYSIS

Shapiro-Wilks test for Normality of Residuals -- alpha-level=0.01

Levenes test for homogeneity of variance(absolute residuals) -- alpha-level=0.05 Use parametric analyses if neither test rejected, otherwise non-parametric analyses.

Shapiro-Wilks	Shapiro-Wilks	Levenes	Levenes	Conclusion
Test Stat	P-value	Test Stat	P-value	
0.928	0.143	0.624	0.652	USE PARAMETRIC TESTS

\*

BASIC ST	JMMARY	STATIS	TICS				
Level	N	Mean	StdDev	StdErr	Coef of Var	95% Conf.I	interval
Ctrl	4	11.94	1.65	0.83	13.83	9.31,	14.57
Dose1	4	12.92	2.29	1.15	17.74	9.27,	16.57
Dose2	4	11.21	2.58	1.29	22.98	7.11,	15.31
Dose3	4	11.81	3.44	1.72	29.16	6.33,	17.29
Dose4	4	13.32	2.02	1.01	15.16	10.10,	16.53
Level		Median	Min	Max	<pre>%of Control(means)</pre>	) %Reduct	ion(means)
Ctrl		11.59	10.48	14.12			
Dose1		13.27	10.12	15.04	108.21	-8.2	1
Dose2		11.81	7.89	13.34	93.89	6.1	.1
Dose3		12.77	6.88	14.81	98.87	1.1	.3
Dose4		13.71	10.52	15.33	111.49	-11.4	.9

\*

PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests

Analysis of Variance (ANOVA) - overall F-test

F-stat P-value 0.48 0.748 Numerator df Denominator df F-stat 4 15

Page 68 of 86

Dunnett - testing each trt mean signif. different than control Williams - test assumes dose-response relationship, testing negative trend Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC

Level	Mean	Dunnett	Isotonic	Williams	Tukey p-values					
		p-value	mean	p-value	Dose1	Dose2	Dose3	Dose4	Dose5	
Ctrl	11.94		12.43			•		•		
Dose1	12.92	0.945	12.43	0.696						
Dose2	11.21	0.980	12.11	0.659	0.861		•			
Dose3	11.81	1.000	12.11	0.678	0.966	0.997	•	•		
Dose4	13.32	0.847	12.11	0.690	0.999	0.750	0.906			

\*\*\*\*\*\*\*\*\*\*\*\*\*\*\*\*\*

NON-PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests

Kruskal-Wallis test - equality among treatment groups

Degrees of Freedom TestStat P-value

2.60 0.627

MannWhit - testing each trt median signif. different from control Jonckheere - test assumes dose-response relationship, testing negative trend

Level	Median	MannWhit p-value	Jonckheere p-value
Ctrl	11.59		
Dose1	13.27	0.678	0.718
Dose2	11.81	0.678	0.330
Dose3	12.77	0.678	0.463
Dose4	13.71	0.494	0.787

DECREASING TREND TEST SUMMARY LOWEST CONCENTRATION SIGNIF. LESS THAN CONTROL Williams >highest dose (no sign. differences) >highest dose (no sign. differences)

\*

PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
Analysis of Variance (ANOVA) - overall F-test

Numerator df Denominator df F-stat P-value

4 15 0.48 0.748

Dunnett - testing each trt mean signif. different than control Williams - test assumes dose-response relationship, testing INCREASING trend Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC

Level	Mean	Dunnett	Isotonic	Williams	Tukey p-values					
		p-value	mean	p-value	Dose1	Dose2	Dose3	Dose4	Dose5	
Ctrl	-11.94		-11.94							
Dose1	-12.92	0.945	-11.98	0.574						
Dose2	-11.21	0.980	-11.98	0.608	0.861	•	•	•		
Dose3	-11.81	1.000	-11.98	0.627	0.966	0.997				
Dose4	-13.32	0.847	-13.32	0.301	0.999	0.750	0.906			

Page 69 of 86

```
NON-PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
    Kruskal-Wallis test - equality among treatment groups
     Degrees of Freedom TestStat P-value
                                             0.627
             4
                                2.60
MannWhit - testing each trt median signif. different from control
Jonckheere - test assumes dose-response relationship, testing INCREASING trend
          Median
                              MannWhit p-value
                                                               Jonckheere p-value
  Ctrl
          -11.59
  Dose1
           -13.27
                                          0.678
                                                                      0.282
            -11.81
                                          0.678
                                                                      0.670
  Dose2
  Dose3
            -12.77
                                          0.678
                                                                       0.537
  Dose4
            -13.71
                                          0.494
                                                                       0.213
 INCREASING TREND TEST SUMMARY LOWEST CONCENTRATION SIGNIF. GREATER THAN
CONTROL
   Williams
                                                      >highest dose (no sign. differences)
   Jonckheere
                                                      >highest dose (no sign. differences)
test for fish screen study - TEST DATA 2 4 D
ANALYSIS RESULTS FOR VARIABLE VAR08 ( M GSI )
TESTS OF ASSUMPTIONS FOR PARAMETRIC ANALYSIS
Shapiro-Wilks test for Normality of Residuals -- alpha-level=0.01
Levenes test for homogeneity of variance(absolute residuals) -- alpha-level=0.05
Use parametric analyses if neither test rejected, otherwise non-parametric
analyses.
  Shapiro-Wilks Shapiro-Wilks Levenes Levenes Conclusion
    Test Stat P-value Test Stat P-value
       0.944
                     0.289
                                        2.115
                                                     0.129
                                                              USE PARAMETRIC TESTS
************************
BASIC SUMMARY STATISTICS
 Level N Mean StdDev StdErr Coef of Var
                                                                    95% Conf.Interval
                                      0.07 11.96

0.11 18.26

0.05 8.43

0.11 19.11

0.05 7.61
  Ctrl 4 1.20 0.14
                                                                     0.97, 1.43
 Dosel 4 1.23 0.22
Dosel 4 1.19 0.10
Dosel 4 1.19 0.23
Dosel 4 1.30 0.10
                                                                       0.87,
                                                                                  1.59
                                                                                  1.35
                                                                        1.03,
                                                                                  1.55
                                                                        0.83,
                                                                         1.14,

        Level
        Median
        Min
        Max
        % of Control(means)
        % Reduction(means)

        Ctrl
        1.25
        1.00
        1.31
        .
        .

        Dosel
        1.30
        0.92
        1.41
        102.29
        -2.29

        Dose2
        1.21
        1.05
        1.29
        98.75
        1.25

        Dose3
        1.21
        0.94
        1.39
        98.75
        1.25

        Dose4
        1.31
        1.19
        1.40
        108.32
        -8.32

*************************
PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
    Analysis of Variance (ANOVA) - overall F-test
     Numerator df Denominator df F-stat
                                                            P-value
                         15
                                              0.33
                                                            0.856
```

Page 70 of 86

Dunnett - testing each trt mean signif. different than control Williams - test assumes dose-response relationship, testing negative trend Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC

Level	Mean	Dunnett	Isotonic	Williams			Tukey p-	values	
		p-value	mean	p-value	Dose1	Dose2	Dose3	Dose4	Dose5
Q1 1	1 00		1 00						
Ctrl	1.20	•	1.22	•	•	•	•	•	•
Dose1	1.23	0.998	1.22	0.651			•	•	
Dose2	1.19	1.000	1.22	0.686	0.996	•	•	•	
Dose3	1.19	1.000	1.22	0.705	0.996	1.000	•	•	
Dose4	1.30	0.817	1.22	0.717	0.972	0.867	0.867		

\*

NON-PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests

Kruskal-Wallis test - equality among treatment groups

Degrees of Freedom TestStat P-value

1.55 0.817

MannWhit - testing each trt median signif. different from control Jonckheere - test assumes dose-response relationship, testing negative trend

Level	Median	MannWhit p-value	Jonckheere p-value
Ctrl	1.25		
Dose1	1.30	0.678	0.718
Dose2	1.21	0.780	0.384
Dose3	1.21	0.889	0.444
Dose4	1.31	0.494	0.703

DECREASING TREND TEST SUMMARY

Williams

Jonckheere

LOWEST CONCENTRATION SIGNIF. LESS THAN CONTROL

>highest dose (no sign. differences)

>highest dose (no sign. differences)

\*

PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
Analysis of Variance (ANOVA) - overall F-test
Numerator df Denominator df F-stat P-value
4 15 0.33 0.856

Dunnett - testing each trt mean signif. different than control Williams - test assumes dose-response relationship, testing INCREASING trend Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC

Level	Mean	Dunnett	Isotonic	Williams			Tukey p-v	values	
		p-value	mean	p-value	Dosel	Dose2	Dose3	Dose4	Dose5
Ctrl	-1.20		-1.20						
Dose1	-1.23	0.998	-1.20	0.586					
Dose2	-1.19	1.000	-1.20	0.620	0.996				
Dose3	-1.19	1.000	-1.20	0.639	0.996	1.000	•		
Dose4	-1.30	0.817	-1.30	0.281	0.972	0.867	0.867		

Page 71 of 86

```
*************************
NON-PARAMETRIC ANALYSES
                                - use alpha-level=0.05 for all tests
     Kruskal-Wallis test - equality among treatment groups
      Degrees of Freedom TestStat P-value
                                              0.817
                                   1.55
MannWhit - testing each trt median signif. different from control
Jonckheere - test assumes dose-response relationship, testing INCREASING trend
          Median
 Level
                               MannWhit p-value
                                                                  Jonckheere p-value
  Ctrl
            -1.25
  Dose1
              -1.30
                                            0.678
                                                                           0.282
  Dose2 -1.21
                                            0.780
                                                                          0.616
  Dose3
             -1.21
                                            0.889
                                                                          0.556
                                            0.494
  Dose4
             -1.31
                                                                          0.297
 INCREASING TREND TEST SUMMARY LOWEST CONCENTRATION SIGNIF. GREATER THAN
CONTROL
   Williams
                                                         >highest dose (no sign. differences)
   Jonckheere
                                                         >highest dose (no sign. differences)
test for fish screen study - TEST DATA 2 4 D
ANALYSIS RESULTS FOR VARIABLE VAR09 ( F tubercle score (median) )
TESTS OF ASSUMPTIONS FOR PARAMETRIC ANALYSIS
Shapiro-Wilks test for Normality of Residuals -- alpha-level=0.01
Levenes test for homogeneity of variance(absolute residuals) -- alpha-level=0.05
Use parametric analyses if neither test rejected, otherwise non-parametric
  Shapiro-Wilks Shapiro-Wilks Levenes Levenes Conclusion
     Test Stat P-value Test Stat P-value
                                                       . NO DATA FOR TEST
**************************
BASIC SUMMARY STATISTICS
Level N Mean StdDev StdErr Coef of Var 95% Conf.Interval

      Ctrl 4
      0.00
      0.00
      0.00
      .

      Dosel 4
      0.00
      0.00
      0.00
      .

      Dose2 4
      0.00
      0.00
      0.00
      .

      Dose3 4
      0.00
      0.00
      0.00
      .

      Dose4 4
      0.00
      0.00
      0.00
      .

        Level
        Median
        Min
        Max
        % of Control(means)
        %Reduction(means)

        Ctrl
        0.00
        0.00
        .
        .

        Dosel
        0.00
        0.00
        .
        .

        Dose2
        0.00
        0.00
        0.00
        .
        .

        Dose3
        0.00
        0.00
        0.00
        .
        .

        Dose4
        0.00
        0.00
        0.00
        .
        .

************************
PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
    Analysis of Variance (ANOVA) - overall F-test
     Numerator df Denominator df F-stat
```

Page 72 of 86

						EPA	MRID Numl	ber 483170	01
		. 1							_
Williams	s - test	assumes	dose-respo	gnif. diff onse relati ble compar	onship,	testing	negative		
Level	Mean	Dunnett p-value		Williams p-value	Dose1		Tukey p- Dose3	values Dose4	Dose5
Ctrl	0.00								
Dosel	0.00	•	•	•	•	•	•	•	•
Dose1 Dose2	0.00	•	•	•	•	•	•	•	•
		•	•	•	•	•	•	•	•
Dose3	0.00	•	•	•	•	•	•	•	•
Dose4	0.00	•	•	•	•	•	•	•	•
*****	*****	*****	******	*****	*****	*****	*****	*****	***
NON-PARA	AMETRIC	ANALYSES	- use a	alpha-level	=0.05 fo	r all te	sts		
				among tre					
		Freedom	TestStat	_	_	τοαρο			
Deč	_	rreedom							
	4		0.00	1.00	U				
N/ T-7] +					. e e	£			
				n signif. d				_	
Jonckhee	ere - te	est assume	es dose-res	sponse rela	tionship	, testin	g negati	ve trend	
Level	Medi	an	MannWhi	t p-value		Jonakh	eere p-v	alue	
Ctrl		00	Mailliwill	ic p varue		OOHCKII	CCIC P V	aruc	
		00		1.000			•		
Dose1							•		
Dose2		00		1.000			•		
Dose3	0.	00		1.000			•		
Dose4	0.	00		1.000			•		
		ND TEST S	SUMMARY	LOWEST C	ONCENTRA'	TION SIG	NIF. LES	S THAN C	ONTROL
Willi					Dose1				
Joncl	kheere				Dose1				
******	******			******			*****	****	***
			_	a-level=0.0		1 tests			
	iysis oi merator			overall F lf F-sta		P-value			
•			•	•	·	·		·	·
Dunnett	- testi	.ng each t	rt mean si	gnif. diff	erent th	an contr	ol		
Williams	s - test	assumes	dose-respo	nse relati	onship,	testing	INCREASI	NG trend	
				ble compar					
-			_	-	•				
Level	Mean	Dunnett	Isotonic	Williams			Tukey p-	values	
		p-value	mean	p-value	Dose1		Dose3		Dose5
			30.22	1 0.1 0.0					
Ctrl	0.00						_		
Dosel	0.00	•	•	•	•	-	•	•	•
Dose1 Dose2	0.00	•	•	•	•	•	•	•	•
DOSEZ	0.00	•	•	•	•	•	•	•	•

Page 73 of 86

						EPA	ואואוט ואנווו	Del 40317	001
Dose3	0.	00 .			•				
Dose4		00 .							
++++++	++++	++++++	*****	****	******	+++++	++++++	++++++	++++
					vel=0.05 for			****	***
				-	treatment gro		515		
			dom TestSt			ups			
Deg	4		0.0		.000				
	7		0.0	J	.000				
					. different f elationship,			SING tre	end
Level		ledian	Mann	Whit p-val	ue	Jonckh	eere p-v	alue	
Ctrl		0.00					•		
Dose1		0.00		1.000			•		
Dose2		0.00		1.000			•		
Dose3		0.00		1.000			•		
Dose4		0.00		1.000			•		
INCREAS CONTROL	ING	TREND TE	ST SUMMARY	LOWES	T CONCENTRATI	ON SIG	NIF. GRE	ATER THA	AN
Willi	ams				Dose1				
Jonck		·e			Dose1				
test for	fis	h screen	study - TES	ST DATA 2	4 D				
ANALYSIS	RES	ULTS FOR	. VARIABLE V	AR10 ( M	tubercle scor	e (med:	ian))		
TESTS OF	' ASS	SUMPTIONS	FOR PARAMET	TRIC ANALY	SIS				
Shapiro-	Wilk	s test f	or Normality	y of Resid	uals alpha	-level:	=0.01		
Levenes	test	for hom	ogeneity of	variance(	absolute resi	duals)	alph	a-level=	=0.05
		ic analy	ses if neith	ner test r	ejected, othe	rwise 1	non-para	metric	
analyses									
_			_		Levenes	Concl	usion		
			-value						
0.	963		0.608	0.683	0.615	USE P	ARAMETRI	C TESTS	
******	****	*****	*****	******	*****	****	******	*****	****
BASIC SU									
_	N	Mean	_	C+dExx	Coef of Va	ac 0.1	E% Conf	Tnton.	1
			StdDev	2.01		L 9:			L
	4	32.00	4.02		12.56		25.60,		
Dose1		31.25	3.88	1.94	12.43		25.07,	37.43	
Dose2		33.63	4.96	2.48	14.74		25.74, 29.98,	41.51	
Dose3		33.13	1.97	0.99	5.96		29.98,	36.27	
Dose4	4	34.50	3.08	1.54	8.93		29.60,	39.40	
Level		Median	Min	Max	%of Control(m	leans)	%Reduc	tion(mea	ang)
Ctrl		31.25	28.00	37.50		,	11.0440	\	,
Dose1		31.00	27.00	36.00	97.66		· 2	34	
Dose1		34.50	27.00	38.50	105.08		-5.		
Dose3		32.50	31.50	36.00	103.52		-3.		
Dose3		34.25	31.00	38.50	103.32		-3. -7.		
POSC4		24.43	31.00	30.30	107.01		- / .	0.1	
*****	***	*****	****	*****	******	****	*****	*****	****

Page 74 of 86

PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
Analysis of Variance (ANOVA) - overall F-test
Numerator df Denominator df F-stat P-value
4 15 0.48 0.749

Dunnett - testing each trt mean signif. different than control Williams - test assumes dose-response relationship, testing negative trend Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC

Level	Mean	Dunnett	Isotonic	Williams			Tukey p-	<i>r</i> alues	
		p-value	mean	p-value	Dose1	Dose2	Dose3	Dose4	Dose5
~. I	20.00		20.00						
Ctrl	32.00	•	32.90	•	•	•		•	
Dose1	31.25	0.995	32.90	0.719			•		
Dose2	33.63	0.925	32.90	0.754	0.891		•		
Dose3	33.13	0.978	32.90	0.772	0.950	1.000	•		
Dose4	34.50	0.748	32.90	0.783	0.732	0.997	0.984		

\*

NON-PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests

Kruskal-Wallis test - equality among treatment groups

Degrees of Freedom TestStat P-value

2.41 0.661

MannWhit - testing each trt median signif. different from control Jonckheere - test assumes dose-response relationship, testing negative trend

Level	Median	MannWhit p-value	Jonckheere p-value
Ctrl	31.25		
Dose1	31.00	0.889	0.386
Dose2	34.50	0.678	0.722
Dose3	32.50	0.491	0.745
Dose4	34.25	0.346	0.879

DECREASING TREND TEST SUMMARY

Williams

Jonckheere

LOWEST CONCENTRATION SIGNIF. LESS THAN CONTROL

>highest dose (no sign. differences)

>highest dose (no sign. differences)

\*

PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
Analysis of Variance (ANOVA) - overall F-test
Numerator df Denominator df F-stat P-value
4 15 0.48 0.749

Dunnett - testing each trt mean signif. different than control Williams - test assumes dose-response relationship, testing INCREASING trend Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC

Level	Mean	Dunnett	Isotonic	Williams			Tukey p-	values	
		p-value	mean	p-value	Dose1	Dose2	Dose3	Dose4	Dose5
Ctrl	-32.00		-31.63						
Dose1	-31.25	0.995	-31.63	0.642					

Page 75 of 86

Dose2 -33	3.63 0	.925	-33.38	0.388	0.891				
Dose3 -33	3.13 0	.978	-33.38	0.401	0.950	1.000			
Dose4 -34			-34.50	0.241	0.732	0.997	0.984		
******								*****	****
NON-PARAMET							ests		
					reatment g	coups			
Degree	s of Fr	eedom 5	TestStat						
	4		2.41	0.	661				
MannWhit -	testing	each tr	t median	signif.	different	from co	ntrol		
Jonckheere	_			_				SING tre	nd
	Median		MannWhi	t p-valu	ıe	Jonckh	neere p-va	alue	
	-31.25								
Dose1				0.889			0.614		
	-34.50			0.678			0.278		
	-32.50			0.491			0.255		
Dose4	-34.25			0.346			0.121		
INCREASING CONTROL	TREND	TEST SUMI	MARY	LOWEST	CONCENTRA	rion sid	SNIF. GRE	ATER THA	N
Williams					shiahe	et doga	(no sign	differ	ancac
	)								
	200				>hiaha	7+ dogo	/no gian		
Jonckhee	ere				>highe	st dose	(no sign	. diller	ences)
Jonckhee		en studv	- TEST	DATA 2 4		st dose	(no sign	. diller	ences)
Jonckhee test for fi	sh scre				ł D	st dose	(no sign	. differ	ences )
Jonckhee	sh scre				ł D	st dose	(no sign	. differ	ences)
Jonckhee test for fi	sh scre	OR VARIA	BLE VAR1	1 ( fec	D Dundity )	st dose	(no sign	. diller	ences )
Jonckhee test for fi ANALYSIS RE TESTS OF AS Shapiro-Wil	sh scre SULTS F SUMPTIO	OR VARIA NS FOR PA	BLE VAR1	l (fec C ANALYS f Residu	l D cundity ) SIS als alpl	na-level	_=0.01		
Jonckhee test for fi ANALYSIS RE TESTS OF AS Shapiro-Wil Levenes tes	sh scre SULTS F SUMPTIO	OR VARIA NS FOR PA for Norm omogenei	BLE VAR1: ARAMETRIC mality of ty of va:	l (fec C ANALYS f Residu riance(a	l D cundity ) SIS als alpl absolute re:	na-level siduals)	=0.01 alph	a-level=	
Jonckhee  test for fi ANALYSIS RE  TESTS OF AS Shapiro-Wil Levenes tes Use paramet	sh scre SULTS F SUMPTIO	OR VARIA NS FOR PA for Norm omogenei	BLE VAR1: ARAMETRIC mality of ty of va:	l (fec C ANALYS f Residu riance(a	l D cundity ) SIS als alpl absolute re:	na-level siduals)	=0.01 alph	a-level=	
Jonckhee test for fi ANALYSIS RETESTS OF AS Shapiro-Wil Levenes test Use paramet analyses.	sh scre SSUMPTIO ks test t for h	OR VARIA NS FOR PA for Normomogenei omogenei lyses if	BLE VAR1  ARAMETRIO  mality of  ty of va:  neither	l (fec C ANALYS f Residu riance(a test re	Deundity ) SIS Hals alphabsolute resisted, other	na-level siduals) nerwise	=0.01 alphanon-param	a-level=	
Jonckhee  test for fi ANALYSIS RE  TESTS OF AS Shapiro-Wil Levenes tes Use paramet analyses. Shapiro-W	sh scre SULTS F SUMPTIO ks test t for h cric ana	OR VARIANT OR VARIANT OR FOR Normogenei lyses if	ARAMETRIO mality of ty of variable neither	1 ( fec C ANALYS f Residu riance(a test re Levenes	Eundity )  SIS  als alpl  absolute reserved, other	na-level siduals) nerwise Concl	=0.01 alph	a-level=	
Jonckhee  test for fi ANALYSIS RE  TESTS OF AS Shapiro-Wil Levenes tes Use paramet analyses. Shapiro-W Test St	sh scre SULTS F SUMPTIO ks test t for h cric ana	OR VARIANT OR VARIANT OR FOR PARTIES OF PARTIES OF THE PROPERTY OF THE PROPERT	ARAMETRIO mality of ty of variable neither	1 ( fec C ANALYS f Residu riance(a test re Levenes est Stat	Levenes Deundity )  Levenes Deundity )	na-level siduals) nerwise Concl	=0.01 alph non-para usion	a-level= metric	
Jonckhee  test for fi ANALYSIS RE  TESTS OF AS Shapiro-Wil Levenes tes Use paramet analyses. Shapiro-W	sh scre SULTS F SUMPTIO ks test t for h cric ana	OR VARIANT OR VARIANT OR FOR Normogenei lyses if	ARAMETRIO mality of ty of variable neither	1 ( fec C ANALYS f Residu riance(a test re Levenes est Stat	Eundity )  SIS  als alpl  absolute reserved, other	na-level siduals) nerwise Concl	=0.01 alphanon-param	a-level= metric	
Jonckhee  test for fi ANALYSIS RE  TESTS OF AS Shapiro-Wil Levenes tes Use paramet analyses. Shapiro-W Test St	sh scre SULTS F SUMPTIO ks test t for h cric ana lilks S	OR VARIANT OR VARIANT OR FOR PARTIES OF THE PROPERTY OF THE PR	ARAMETRIO mality o ty of va: neither ilks	1 ( fec C ANALYS f Residu riance(a test re Levenes est Stat 0.391	Levenes P-value 0.812	na-level siduals) nerwise Concl USE I	=0.01 alphanon-parameters 	a-level= metric C TESTS	0.05
Jonckhee  test for fi ANALYSIS RE  TESTS OF AS Shapiro-Wil Levenes tes Use paramet analyses. Shapiro-W Test St 0.964	sh scre SULTS F SUMPTIO ks test t for h cric ana Tilks S at	OR VARIANT	ARAMETRIO mality o ty of va: neither ilks	1 ( fec C ANALYS f Residu riance(a test re Levenes est Stat 0.391	Levenes P-value 0.812	na-level siduals) nerwise Concl USE I	=0.01 alphanon-parameters 	a-level= metric C TESTS	0.05
Jonckhee  test for fi ANALYSIS RE  TESTS OF AS Shapiro-Wil Levenes tes Use paramet analyses. Shapiro-W Test St 0.964  ***********************************	sh scre SULTS F SUMPTIO ks test t for h cric ana lilks S at t xxxxxxxxxxxxxxxxxxxxxxxxxxxxxxxxx	OR VARIANT	ARAMETRIC mality of ty of variable neither ilks	1 ( fee C ANALYS f Residu riance(a test re Levenes est Stat 0.391	Levenes P-value 0.812	na-level siduals) nerwise Concl USE E	=0.01 alphonon-parametrical constants on Parametrical constants on the constant of the cons	a-level= metric C TESTS *****	0.05
Jonckhee  test for fi ANALYSIS RE  TESTS OF AS Shapiro-Wil Levenes tes Use paramet analyses. Shapiro-W Test St 0.964  ***********************************	sh scre SULTS F SUMPTIO ks test t for h cric ana lilks S at  ******* MY STAT Mean	OR VARIANT	ARAMETRIC mality of ty of variation in the residue of the residue	1 ( fee C ANALYS f Residuriance(a test re Levenes est Stat 0.391 ******	Levenes P-value 0.812	na-level siduals) nerwise Concl USE E	=0.01 alphonon-parametric PARAMETRIC ********	a-level= metric C TESTS ******	0.05
Jonckhee  test for fi ANALYSIS RE  TESTS OF AS Shapiro-Wil Levenes tes Use paramet analyses. Shapiro-W Test St 0.964  ***********************************	sh scre SULTS F SUMPTIO ks test t for h cric ana lilks S at  *******  RY STAT Mean 29.8	OR VARIANT	ARAMETRIC mality of ty of value neither liks To see the see th	1 ( fee C ANALYS f Residuriance(a test re Levenes est Stat 0.391 *******	Levenes P-value 0.812  Coef of 121.47	na-level siduals) nerwise Concl USE E	=0.01 alphanon-paramusion PARAMETRIC ********  25% Conf.: 19.62,	a-level= metric  C TESTS  ******  Interval 39.98	0.05
Jonckhee  test for fi ANALYSIS RE  TESTS OF AS Shapiro-Wil Levenes tes Use paramet analyses. Shapiro-W Test St 0.964  *********** BASIC SUMMA Level N Ctrl 4 Dosel 4	sh scre SULTS F SUMPTIO ks test t for h cric ana lilks S at  ******* RY STAT Mean 29.8 29.0	OR VARIANT	ARAMETRIC mality of ty of variation of variation of the meither liks To the meither with the meither t	C ANALYS f Residuriance(a test re Levenes est Stat 0.391 ******* StdErr 3.20 3.23	Levenes P-value 0.812  Coef of v 21.47 22.28	na-level siduals) nerwise Concl USE E	=0.01 alphonon-parametrical control contro	a-level= metric  C TESTS  ******  Interval 39.98 39.31	0.05
Jonckhee  test for fi ANALYSIS RE  TESTS OF AS Shapiro-Wil Levenes tes Use paramet analyses. Shapiro-W Test St 0.964  *********** BASIC SUMMA Level N Ctrl 4 Dosel 4 Dose2 4	sh scre SULTS F SUMPTIO ks test t for h cric ana lilks S at  ****** RY STAT Mean 29.8 29.0 25.9	OR VARIANT	ARAMETRIC mality of ty of variather ilks  ********  ev 40 47 27	1 ( fee C ANALYS f Residuriance(a test re Levenes est Stat 0.391 ******** StdErr 3.20 3.23 2.64	Levenes P-value 0.812  Coef of v 21.47 22.28 20.34	na-level siduals) nerwise Concl USE E	=0.01 alphonon-paramusion  PARAMETRIC  ********  95% Conf.: 19.62, 18.74, 17.53,	a-level= metric  C TESTS  ******  Interval 39.98 39.31 34.32	0.05
Jonckhee  test for fi ANALYSIS RE  TESTS OF AS Shapiro-Wil Levenes tes Use paramet analyses. Shapiro-W Test St 0.964  *********** BASIC SUMMA Level N Ctrl 4 Dosel 4 Dose2 4 Dose3 4	sh scre SSUMPTIO ks test t for h cric ana lilks S at  ****** *** **** Mean 29.8 29.0 25.9 22.8	OR VARIANT	ARAMETRIC mality of ty of variation of variation of variation of the results.  ***********  ev 40 47 27 72	1 (fee C ANALYS f Residuriance(a test re Levenes est Stat 0.391 ******** StdErr 3.20 3.23 2.64 4.36	Levenes P-value 0.812  Coef of v 21.47 22.28	na-level siduals) nerwise Concl USE E	=0.01 alphonon-parametrical control contro	a-level= metric C TESTS ******* Interval 39.98 39.31 34.32 36.76	0.05
Jonckhee  test for fi ANALYSIS RE  TESTS OF AS Shapiro-Wil Levenes tes Use paramet analyses. Shapiro-W Test St 0.964  *********** BASIC SUMMA Level N Ctrl 4 Dosel 4 Dose2 4 Dose3 4	sh scre SULTS F SUMPTIO ks test t for h cric ana lilks S at  ****** RY STAT Mean 29.8 29.0 25.9	OR VARIANT	ARAMETRIC mality of ty of variation of variation of variation of the results.  ***********  ev 40 47 27 72	1 ( fee C ANALYS f Residuriance(a test re Levenes est Stat 0.391 ******** StdErr 3.20 3.23 2.64	Levenes P-value 0.812  Coef of v 21.47 22.28 20.34 38.14	na-level siduals) nerwise Concl USE E	=0.01 alphonon-parametrical control contro	a-level= metric  C TESTS  ******  Interval 39.98 39.31 34.32	0.05
Jonckhee  test for fi ANALYSIS RE  TESTS OF AS Shapiro-Wil Levenes tes Use paramet analyses. Shapiro-W Test St 0.964  *********** BASIC SUMMA Level N Ctrl 4 Dosel 4 Dose2 4 Dose3 4	sh scre SSUMPTIO ks test t for h cric ana lilks S at  ****** *** **** Mean 29.8 29.0 25.9 22.8	OR VARIANT	ARAMETRIC mality of ty of variable in the results and results are	1 (fee C ANALYS f Residuriance(a test re Levenes est Stat 0.391 ******** StdErr 3.20 3.23 2.64 4.36 2.98	Levenes P-value 0.812  Coef of v 21.47 22.28 20.34 38.14	na-level siduals) nerwise Concl USE E	=0.01 alphonon-parametrical control contro	a-level= metric C TESTS ******** Interval 39.98 39.31 34.32 36.76 29.20	0.05 ***
Jonckhee  test for fi ANALYSIS RE  TESTS OF AS Shapiro-Wil Levenes tes Use paramet analyses. Shapiro-W Test St 0.964  **********  BASIC SUMMA Level N Ctrl 4 Dosel 4	sh scre SULTS F SUMPTIO ks test t for h ric ana lilks S at  ****** Mean 29.8 29.0 25.9 22.8 19.7	OR VARIANT	ARAMETRIC mality of ty of variable in the results i	1 (fee C ANALYS f Residuriance(a test re Levenes est Stat 0.391 ******** StdErr 3.20 3.23 2.64 4.36 2.98	Levenes P-value 0.812  Coef of V 21.47 22.28 20.34 38.14 30.29	na-level siduals) nerwise Concl USE E	=0.01 alphonon-parametrical control contro	a-level= metric  C TESTS  ******  Interval	0.05 ***
Jonckhee  test for fi ANALYSIS RE  TESTS OF AS Shapiro-Wil Levenes tes Use paramet analyses. Shapiro-W Test St 0.964  **********  BASIC SUMMA Level N Ctrl 4 Dosel 4 Dosel 4 Dose2 4 Dose3 4 Dose4 4  Level	sh scre SULTS F SUMPTIO ks test t for h ric ana Tilks S at  ****** MEAN 29.8 29.0 25.9 22.8 19.7 Media	OR VARIANT	ARAMETRIC mality of ty of variable in the results and results are	1 (fee C ANALYS f Residuriance(a test re Levenes est Stat 0.391 ******** StdErr 3.20 3.23 2.64 4.36 2.98	Levenes P-value 0.812  Coef of 121.47 22.28 20.34 38.14 30.29	na-level siduals) nerwise Concl USE E	=0.01 alphanon-paramusion PARAMETRIC ************************************	a-level= metric  C TESTS  ******  Interval     39.98     39.31     34.32     36.76     29.20  tion(mea	0.05 ***
Jonckhee  test for fi ANALYSIS RE  TESTS OF AS Shapiro-Wil Levenes tes Use paramet analyses. Shapiro-W Test St 0.964  **********  BASIC SUMMA Level N Ctrl 4 Dosel 4	sh scre SULTS F SUMPTIO ks test t for h cric ana Tilks S tat  ****** RY STAT Mean 29.8 29.0 25.9 22.8 19.7 Media 28.3 29.0	OR VARIANT	ARAMETRIC mality of ty of variation in the residue of the residue	1 ( fee C ANALYS f Residuriance(a test re Levenes est Stat 0.391 ******** StdErr 3.20 3.23 2.64 4.36 2.98 Max % 38.60 35.10	Levenes P-value 0.812  Coef of 21.47 22.28 20.34 38.14 30.29  Sof Control . 97.40	na-level siduals) nerwise Concl USE E	=0.01 alphonon-paramusion  PARAMETRIC  ********  95% Conf.: 19.62, 18.74, 17.53, 8.99, 10.20,  %Reduc 2.1	a-level= metric  C TESTS  *******  Interval     39.98     39.31     34.32     36.76     29.20  tion(mea	0.05 ***
Jonckhee  test for fi ANALYSIS RE  TESTS OF AS Shapiro-Wil Levenes tes Use paramet analyses. Shapiro-W Test St 0.964  ********** BASIC SUMMA Level N Ctrl 4 Dosel 4 Dosel 4 Dosel 4 Dosel 4 Dosel 4 Dosel 4 Level Ctrl Ctrl	sh scre SULTS F SUMPTIO ks test t for h ric ana Tilks S at  ****** ARY STAT Mean 29.8 29.0 25.9 22.8 19.7 Media 28.3	OR VARIANT	ARAMETRIC mality of ty of variation of varia	1 ( fee C ANALYS f Residuriance(a test re Levenes est Stat 0.391 ******** StdErr 3.20 3.23 2.64 4.36 2.98 Max %	Levenes P-value 0.812  *********  Coef of 121.47 22.28 20.34 38.14 30.29  Sof Control	na-level siduals) nerwise Concl USE E	=0.01 alphanon-paramusion  PARAMETRIC  ********  95% Conf. 19.62, 18.74, 17.53, 8.99, 10.20,  %Reduc	a-level= metric  C TESTS  ******  Interval     39.98     39.31     34.32     36.76     29.20  tion(means)  60 00	0.05

Page 76 of 86

```
*************************
PARAMETRIC ANALYSES

    use alpha-level=0.05 for all tests

    Analysis of Variance (ANOVA) - overall F-test
    Numerator df Denominator df F-stat
                                                     P-value
                      15
                                        1.61
                                                     0.223
Dunnett - testing each trt mean signif. different than control
Williams - test assumes dose-response relationship, testing negative trend
Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC
                Dunnett Isotonic Williams
Level
       Mean
                                                              Tukey p-values
                p-value mean p-value Dose1 Dose2 Dose3 Dose4

      Ctrl
      29.80
      .
      .
      .
      .
      .

      Dosel
      29.03
      0.999
      29.03
      0.513
      .
      .

      Dose2
      25.93
      0.826
      25.93
      0.271
      0.962
      .

      Dose3
      22.88
      0.418
      22.88
      0.105
      0.693
      0.965

      Dose4
      19.70
      0.146
      19.70
      0.031
      0.322
      0.683

                                              0.322 0.683
                                                                 0.959
************************
NON-PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
    Kruskal-Wallis test - equality among treatment groups
     Degrees of Freedom TestStat P-value
                              5.36
                                        0.252
MannWhit - testing each trt median signif. different from control
Jonckheere - test assumes dose-response relationship, testing negative trend
Level
         Median
                          MannWhit p-value
                                                        Jonckheere p-value
          28.30
 Ctrl
                                     0.780
 Dose1
           29.05
                                                              0.332
           27.15
 Dose2
                                     0.678
                                                              0.210
 Dose3
           22.85
                                     0.346
                                                               0.088
 Dose4
           18.35
                                     0.156
                                                               0.011
 DECREASING TREND TEST SUMMARY LOWEST CONCENTRATION SIGNIF. LESS THAN CONTROL
   Williams
                                                Dose4
   Jonckheere
                                                Dose4
************************
PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
    Analysis of Variance (ANOVA) - overall F-test
    Numerator df Denominator df F-stat
                                                    P-value
                                        1.61
Dunnett - testing each trt mean signif. different than control
Williams - test assumes dose-response relationship, testing INCREASING trend
Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC
Level
                Dunnett Isotonic Williams
                                                              Tukey p-values
                p-value mean p-value Dose1 Dose2 Dose3 Dose4 Dose5
Ctrl -29.80 .
                          -25.47
                                                        .
```

Page 77 of 86

NON-PARAMET Kruskal Degree MannWhit -	2.93 0.8 2.88 0.4 2.70 0.1 2************************************	26 -25.47 18 -25.47 46 -25.47  ******** SES - use est - equali dom TestSt 5.3  ach trt medi	0.907 0.919 0.926 ************************************	0.693 0 0.322 0 ***********************************		
Level Ctrl Dosel Dose2 Dose3	Median		hit p-value 0.780 0.678 0.346 0.156		onckheere p-v 0.668 0.790 0.912 0.989	
CONTROL Williams Jonckhee test for fi ANALYSIS RE	re sh screen	ST SUMMARY  study - TES  VARIABLE VA	T DATA 2 4 R12 ( fert	>highest >highest D	_	ATER THAN . differences) . differences)
Shapiro-Wil Levenes tes	ks test for home ric analystilks Shapat P	or Normality ogeneity of ses if neith piro-Wilks -value	of Residua variance(ak er test re Levenes Test Stat	als alpha- psolute resid jected, other Levenes P-value	level=0.01 uals) alph wise non-para Conclusion USE PARAMETRI	metric
********** BASIC SUMMA Level N Ctrl 4 Dosel 4 Dose2 4 Dose3 4 Dose4 4  Level Ctrl Dose1 Dose2 Dose3 Dose4			StdErr 0.25 1.54 0.77 0.39 1.34	Coef of Var 0.52 3.19 1.60 0.80 2.79 of Control(me 99.36 99.17 100.77 98.84	96.13, 91.41, 93.68, 96.43, 91.54,	Interval 97.72 101.19 98.57 98.92 100.06 tion(means) 64 83 77

Page 78 of 86

Tukey p-values

p-value Dose1 Dose2 Dose3 Dose4

\* PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests Analysis of Variance (ANOVA) - overall F-test Numerator df Denominator df F-stat 4 15 0.56 0.697 Dunnett - testing each trt mean signif. different than control Williams - test assumes dose-response relationship, testing negative trend Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC Level Mean Dunnett Isotonic Williams Tukey p-values p-value mean p-value Dose1 Dose2 Dose3 Dose4 Dose5 

 Ctrl
 96.93
 .
 96.93
 .
 .

 Dosel
 96.30
 0.975
 96.70
 0.515
 .

 Dose2
 96.13
 0.943
 96.70
 0.547
 1.000

 Dose3
 97.68
 0.954
 96.70
 0.565
 0.862

 Dose4
 95.80
 0.839
 95.80
 0.206
 0.206

 0.804 0.999 0.677 \* NON-PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests Kruskal-Wallis test - equality among treatment groups Degrees of Freedom TestStat P-value 3.93 0.415 MannWhit - testing each trt median signif. different from control Jonckheere - test assumes dose-response relationship, testing negative trend Level Median MannWhit p-value Jonckheere p-value Ctrl 97.10 97.75 0.343 0.877 Dose1 Dose2 96.60 0.678 0.385 97.70 0.235 Dose3 0.825 95.30 0.346 0.474 Dose4 DECREASING TREND TEST SUMMARY LOWEST CONCENTRATION SIGNIF. LESS THAN CONTROL Williams >highest dose (no sign. differences) Jonckheere >highest dose (no sign. differences) \* PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests Analysis of Variance (ANOVA) - overall F-test Numerator df Denominator df F-stat 4 0.56 0.697 15 Dunnett - testing each trt mean signif. different than control Williams - test assumes dose-response relationship, testing INCREASING trend Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC Mean Dunnett Isotonic Williams

Page 79 of 86

p-value mean

Level

Version: 22 September 2011

Dose5

Ctrl	-96.93		-96.45	•				•	
Dose1			5 -96.45	0.717	7.			•	
Dose2	-96.13	0.94	3 -96.45		1.000	•		•	
			4 -96.74		0.862	0.804		•	
Dose4	-95.80	0.83	9 -96.74	0.704	0.996	0.999	0.677	•	
*****	****	*****	*****	*****	*****	*****	*****	*****	***
NON-PARA	METRIC Z	ANALYS	ES - use	alpha-lev	rel=0.05 for	r all te	sts		
					reatment gi				
Deg:	rees of	Freed	om TestSt	at P-va	alue				
	4		3.9	3 0.	. 415				
MannWhit.	- test	ing ea	ch trt medi	an signif.	. different	from co	ntrol		
					elationship			ASING tre	nd
Level	Media	an.	MannW	hit p-valı	10	Tonakh	eere p-v	72 ] 110	
Ctrl	-97.1		Mailliw	iic p varc		UOIICKII	eere b (	varue	
Dosel				0.343			0.123		
Dose2	-96.			0.678			0.615		
Dose3	-97.	70		0.235			0.175		
Dose4	-95.X			0.346			0.526		
20201	,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,			0.010			0.020		
	ING TRE	ND TES	T SUMMARY	LOWEST	CONCENTRA!	rion sig	NIF. GRE	EATER THA	N
CONTROL Willi	am a				. 1. 1. 1.		, .	3:55	0000
								n. differ n. differ	
Jonck:								n. differ	
Jonck	heere	creen	studv - TES	T DATA 2 4	>highes				
Jonck test for	heere fish so		study – TES VARIABLE VA		>highes	st dose	(no sigr		
Jonck test for	heere fish so				>highes	st dose	(no sigr		
Jonck test for ANALYSIS	heere fish so RESULTS	S FOR		R13 ( F V	>highes 1 D /TG re-analy	st dose	(no sigr		
Jonck test for ANALYSIS TESTS OF	heere fish so RESULTS	S FOR	VARIABLE VA	R13 (F V	>highes 1 D /TG re-analy	st dose ysis (ng	(no sigr		
Jonck test for ANALYSIS TESTS OF Shapiro-	heere fish so RESULT: ASSUMP: Wilks to	FOR TIONS TO SET	VARIABLE VA FOR PARAMET r Normality	R13 ( F N RIC ANALYS of Residu	>highes I D /TG re-analy SIS uals alph	st dose ysis (ng na-level	(no sign /mL) ) =0.01	n. differ	ence
Jonck Lest for ANALYSIS LESTS OF Shapiro- Levenes	fish so RESULT: ASSUMPT Wilks to test for	FOR TIONS : est for homogeneous contracts to the state of	VARIABLE VA FOR PARAMET r Normality geneity of	R13 (F N RIC ANALYS of Residu variance(a	>highes  I D  TG re-analy  SIS  Lals alphabsolute res	st dose ysis (ng na-level siduals)	(no sign /mL) ) =0.01 alph	n. differ na-level=	ence
Jonck test for ANALYSIS TESTS OF Shapiro- Levenes Jse para	fish so RESULTS  ASSUMPT Wilks to test for metric a	FOR TIONS : est for homogeneous contracts to the state of	VARIABLE VA FOR PARAMET r Normality geneity of	R13 (F N RIC ANALYS of Residu variance(a	>highes I D /TG re-analy SIS uals alph	st dose ysis (ng na-level siduals)	(no sign /mL) ) =0.01 alph	n. differ na-level=	ence
Jonck Lest for ANALYSIS LESTS OF Chapiro- Levenes Use paramanalyses	fish so RESULTS  ASSUMPS Wilks to test formetric a	FOR TIONS : est for homoganalys	VARIABLE VA FOR PARAMET r Normality geneity of es if neith	R13 (F N RIC ANALYS of Residu variance(a er test re	>highes  I D  TG re-analy  SIS  Lals alphabsolute res	ysis (ng na-level siduals) nerwise	(no sign /mL) ) =0.01 alph non-para	n. differ na-level=	ence
Jonck Lest for ANALYSIS LESTS OF Chapiro Levenes Jse parai	fish so RESULTS  ASSUMPS Wilks to test formetric a	FOR TIONS : est for homoganalys	VARIABLE VA FOR PARAMET r Normality geneity of es if neith iro-Wilks	R13 (F N RIC ANALYS of Residu variance(a er test re	>highes  I D  TG re-analy  SIS  Lals alph  Absolute resepted, oth  Levenes	ysis (ng na-level siduals) nerwise	(no sign /mL) ) =0.01 alph non-para	n. differ na-level=	ence
Jonck Lest for LNALYSIS CESTS OF Chapiro Levenes Jse parai Lnalyses Shapir Test	fish so RESULTS  ASSUMPT Wilks to test for metric a co-Wilks Stat	FIONS: PIONS: Poor homographic Shap Poor poor poor poor poor poor poor poor	VARIABLE VA FOR PARAMET r Normality geneity of es if neith iro-Wilks value	R13 (F N RIC ANALYS of Residu variance(a er test re Levenes Test Stat	>highes  I D  TG re-analy  SIS  Lals alph  Absolute res  Ejected, oth	ysis (ng na-level siduals) nerwise	(no sign /mL) ) =0.01 alph non-para	n. differ na-level= ametric	ence
Jonck Lest for ANALYSIS TESTS OF Chapiro- Levenes Jse paramanalyses Shapiro Test 0.	fish so RESULTS  ASSUMPT Wilks to test for metric a co-Wilks Stat 958	FIONS: est for homographic homographic homographic shap p-0	VARIABLE VA FOR PARAMET r Normality geneity of es if neith iro-Wilks value .511	RI3 (F N RIC ANALYS of Residu variance(a er test re Levenes Test Stat 0.324	>highes  I D  TG re-analy  SIS  Lals alph  Absolute reserved, oth  Levenes  P-value  0.858	ysis (ng na-level siduals) nerwise n Concl	(no sign  /mL) )  =0.01  alph non-para  usion  ARAMETRI	na-level= ametric	ence
Jonck Test for ANALYSIS TESTS OF Shapiro- Levenes Jse paral analyses Shapir Test 0.	fish so RESULT:  ASSUMPT Wilks to test formetric a .  O-Wilks Stat 958	FIONS: PEST FOR TIONS: PEST FOR THE PEST FOR	VARIABLE VA FOR PARAMET r Normality geneity of es if neith iro-Wilks value .511 *********	RI3 (F N RIC ANALYS of Residu variance(a er test re Levenes Test Stat 0.324	>highes  I D  TTG re-analy  SIS  Lals alph  Absolute reserved, oth  Levenes  D-value	ysis (ng na-level siduals) nerwise n Concl	(no sign  /mL) )  =0.01  alph non-para  usion  ARAMETRI	na-level= ametric	ence
Jonck Lest for ANALYSIS TESTS OF Shapiro- Levenes Jse paral analyses Shapir Test 0.	fish so RESULTS  ASSUMPT Wilks to test for metric a co-wilks Stat 958  **********************************	FIONS: PEST FOR TIONS: PEST FOR TIONS: Shap PEST FOR TIONS TATIST	VARIABLE VA FOR PARAMET r Normality geneity of es if neith iro-Wilks value .511 **********************************	RI3 (F N RIC ANALYS of Residu variance(a er test re Levenes Test Stat 0.324 ***********	>highes  I D  ITG re-analy  ISIS  Lals alphabsolute reserved, oth  Levenes  P-value  0.858	ysis (ng na-level siduals) nerwise Concli USE P	(no sign  /mL) )  =0.01  alph non-para  usion  ARAMETRI  *******	na-level= ametric  IC TESTS  *******	0.05
Jonck Lest for ANALYSIS TESTS OF Shapiro- Levenes Jse paral analyses Shapir Test 0. ******** BASIC SUI Level	fish so RESULTS  ASSUMPT Wilks to test for metric a .  o-Wilks Stat 958  ***********************************	FIONS: PEST for homograph period Shap Period TATIST	VARIABLE VA FOR PARAMET r Normality geneity of es if neith iro-Wilks value .511 ************* ICS StdDev	RI3 (F N RIC ANALYS of Residu variance(a er test re Levenes Test Stat 0.324 *********** StdErr	>highes  I D  TG re-analy  SIS  Lals alphabsolute reserved, oth  Levenes  P-value  0.858	ysis (ng na-level siduals) nerwise Concli USE P	(no sign  /mL) )  =0.01  alph non-para  usion  ARAMETRI  *******	na-level= ametric  IC TESTS  ********	0.05
Jonck Lest for ANALYSIS TESTS OF Shapiro- Levenes Jse paramalyses Shapir Test 0. ******** BASIC SUI Level Ctrl	fish so RESULTS  ASSUMPT Wilks to test for metric a . o-Wilks Stat 958  ***********************************	FIONS: PEST FOR TIONS: PEST FOR THE PEST FOR	VARIABLE VA FOR PARAMET r Normality geneity of es if neith iro-Wilks value .511 ************ ICS StdDev 26.56	RI3 (F N RIC ANALYS of Residu variance(a er test re Levenes Test Stat 0.324 *********** StdErr 13.28	>highes  I D  /TG re-analy  SIS  Lals alphabsolute reserved, oth  Levenes  P-value  0.858	ysis (ng na-level siduals) nerwise Concli USE PA	(no sign  /mL) )  =0.01  alph non-para  usion  ARAMETRI  ******  5% Conf. 17.53,	na-level= ametric  IC TESTS  ********  Interval 102.06	0.05
Jonck: Lest for ANALYSIS LESTS OF Chapiro Levenes Jse paramalyses Shapir Test 0.  ******** BASIC SU Level Ctrl Dosel	fish so RESULTS  ASSUMPT Wilks to test for metric a . o-Wilks Stat 958  ***********************************	FOR TIONS To sest for homograph analys Shap P-0 CONTRACTION TATIST SEAR P.80 4.85	VARIABLE VA FOR PARAMET r Normality geneity of es if neith iro-Wilks value .511  ************ ICS StdDev 26.56 30.98	RI3 (F N RIC ANALYS of Residu variance(a er test re Levenes Test Stat 0.324 ********** StdErr 13.28 15.49	>highes  A D  /TG re-analy  SIS  pals alphabsolute reserved, oth  Levenes  P-value  0.858  **********************************	ysis (ng na-level: siduals) nerwise n Concli USE P.	(no sign  /mL) )  =0.01  alph non-para  usion  ARAMETRI  ******  5% Conf. 17.53, 25.54,	na-level= ametric  IC TESTS  *******  Interval 102.06 124.15	0.05
Jonck Lest for ANALYSIS TESTS OF Shapiro- Levenes Jse paramanalyses Shapir Test 0. ******* BASIC SU Level Ctrl Dose1 Dose2	fish so RESULT:  ASSUMP: Wilks to test for metric a . o-Wilks Stat 958  ******* MMARY STAN M6 4 55 4 74 4 55	FOR TIONS TO SET TO	VARIABLE VA FOR PARAMET r Normality geneity of es if neith iro-Wilks value .511  ********** ICS StdDev 26.56 30.98 18.79	RI3 (F N RIC ANALYS of Residu variance(a er test re Levenes Test Stat 0.324  *********  StdErr 13.28 15.49 9.40	>highes  A D  /TG re-analy  SIS  pals alphabsolute reserved, oth  Levenes P-value 0.858  ***********  Coef of V 44.42 41.40 33.76	ysis (ng na-level: siduals) nerwise n Concl: USE P.	(no sign  /mL) )  =0.01  alph non-para  usion  ARAMETRI  ******  5% Conf. 17.53, 25.54, 25.77,	na-level= ametric  IC TESTS  *******  Interval 102.06 124.15 85.57	0.05
Jonck Test for ANALYSIS TESTS OF Chapiro- Levenes Jse paramanalyses Shapiro Test O. Test O. Test Ctrl Dosel Dose2 Dose3	fish so RESULTS  ASSUMPT Wilks to test for metric a co-Wilks Stat 958  ******* MMARY ST N Me 4 59 4 74 4 59 4 58	FOR TIONS : est for homograph of the control of the	VARIABLE VA FOR PARAMET r Normality geneity of es if neith iro-Wilks value .511  ********** ICS StdDev 26.56 30.98 18.79 33.69	RI3 (F N RIC ANALYS of Residu variance(a er test re Levenes Test Stat 0.324  *********  StdErr 13.28 15.49 9.40 16.84	>highes  A D  /TG re-analy  SIS  pals alphabsolute reserved, oth  Levenes P-value 0.858  **********  Coef of V 44.42 41.40 33.76 57.39	ysis (ng na-level: siduals) nerwise n Concl: USE P.	(no sign  /mL) )  =0.01  alph non-para  usion  ARAMETRI  *******  5% Conf. 17.53, 25.54, 25.77, 5.09,	na-level= ametric  IC TESTS  *******  Interval 102.06 124.15 85.57 112.31	0.05
Jonck Lest for NALYSIS LESTS OF Chapiro Levenes Jse paral Levenes Shapir Test O.  ****** ASIC SU Level Ctrl Dose1 Dose2	fish so RESULTS  ASSUMPT Wilks to test for metric a co-Wilks Stat 958  ******* MMARY ST N Me 4 59 4 74 4 59 4 58	FOR TIONS : est for homograph of the control of the	VARIABLE VA FOR PARAMET r Normality geneity of es if neith iro-Wilks value .511  ********** ICS StdDev 26.56 30.98 18.79	RI3 (F N RIC ANALYS of Residu variance(a er test re Levenes Test Stat 0.324  *********  StdErr 13.28 15.49 9.40	>highes  A D  /TG re-analy  SIS  pals alphabsolute reserved, oth  Levenes P-value 0.858  ***********  Coef of V 44.42 41.40 33.76	ysis (ng na-level: siduals) nerwise n Concl: USE P.	(no sign  /mL) )  =0.01  alph non-para  usion  ARAMETRI  *******  5% Conf. 17.53, 25.54, 25.77, 5.09,	na-level= ametric  IC TESTS  *******  Interval 102.06 124.15 85.57	0.05
Jonck Lest for ANALYSIS LESTS OF Shapiro-Levenes Jse paramalyses Shapiro Test O.  Example 1 Ctrl Dose1 Dose2 Dose3 Dose4	fish so RESULTS  ASSUMPT Wilks to test formetric a co-Wilks Stat 958  ****** MMARY ST N Mo 4 59 4 74 4 59 4 99	S FOR TIONS TO SET TO TO TO THE TO T	VARIABLE VA FOR PARAMET r Normality geneity of es if neith iro-Wilks value .511  ********** ICS StdDev 26.56 30.98 18.79 33.69	RI3 (F N RIC ANALYS of Residu variance(a er test re Levenes Test Stat 0.324  *********  StdErr 13.28 15.49 9.40 16.84 13.82	>highes  A D  /TG re-analy  SIS  pals alphabsolute reserved, oth  Levenes P-value 0.858  **********  Coef of V 44.42 41.40 33.76 57.39	st dose  ysis (ng  na-level: siduals) nerwise :  Concl: USE P. ********  Var 9:	(no sign  /mL) )  =0.01  alph non-para  usion  ARAMETRI  *******  5% Conf. 17.53, 25.54, 25.77, 5.09, 48.51,	na-level= ametric  IC TESTS  *******  Interval 102.06 124.15 85.57 112.31 136.45	ence
Jonck Test for ANALYSIS TESTS OF Chapiro- Levenes Jse paramanalyses Shapiro Test O. Test O. Test Ctrl Dosel Dose2 Dose3	fish so RESULTS  ASSUMPT Wilks to test formetric a co-Wilks Stat 958  ***********************************	S FOR TIONS TO SET TO TO THE T	VARIABLE VA FOR PARAMET r Normality geneity of es if neith iro-Wilks value .511  ************* ICS StdDev 26.56 30.98 18.79 33.69 27.63 Min	RI3 (F N RIC ANALYS of Residu variance(a er test re Levenes Test Stat 0.324  *********  StdErr 13.28 15.49 9.40 16.84 13.82	>highes  A D  /TG re-analy  SIS  pals alphabsolute reserved, oth  Levenes P-value 0.858  **********  Coef of V 44.42 41.40 33.76 57.39 29.88	st dose  ysis (ng  na-level: siduals) nerwise :  Concl: USE P. ********  Var 9:	(no sign  /mL) )  =0.01  alph non-para  usion  ARAMETRI  *******  5% Conf. 17.53, 25.54, 25.77, 5.09, 48.51,	na-level= ametric  IC TESTS  *******  Interval 102.06 124.15 85.57 112.31 136.45	ence
Jonck  Lest for ANALYSIS  LESTS OF Chapiro- Levenes  Jose paramanalyses Shapir  Test  0.  Level Ctrl Dose1 Dose2 Dose3 Dose4  Level Ctrl	fish so RESULTS  ASSUMPT Wilks to test formetric a .  O-Wilks Stat 958  ******* MMARY STATE A .  4 59 4 79 4 59 4 92  Med 60	S FOR TIONS to sest for homograph of the property of the prope	VARIABLE VA FOR PARAMET r Normality geneity of es if neith iro-Wilks value .511  ************** ICS StdDev 26.56 30.98 18.79 33.69 27.63  Min 27.03	R13 (F N RIC ANALYS of Residu variance(a er test re Levenes Test Stat 0.324  *********  StdErr 13.28 15.49 9.40 16.84 13.82  Max 91.27	>highes  A D  /TG re-analy  SIS  Lals alphabsolute reserved, oth  Levenes P-value 0.858  **********  Coef of V 44.42 41.40 33.76 57.39 29.88  A Control	st dose  ysis (ng na-level siduals) nerwise  Concli USE P. *******  Var 9:	(no sign  /mL) )  =0.01 alph non-para  usion  ARAMETRI  *******  5% Conf. 17.53, 25.54, 25.77, 5.09, 48.51, %Reduce	na-level= ametric  IC TESTS  *******  Interval 102.06 124.15 85.57 112.31 136.45  ction(mea	ence
Jonck test for ANALYSIS TESTS OF Shapiro- Levenes Use paramanalyses Shapir Test 0. ******** BASIC SU Level Ctrl Dose1 Dose2 Dose3 Dose4 Level Ctrl Dose1	fish so RESULTS  ASSUMPT Wilks to test formetric a .  O-Wilks Stat 958  ******* MMARY STATE A .  4 59 4 79 4 59 4 92  Med 60	S FOR TIONS to sest for homograph analys Shap P-0 STATIST SEAN P-1 STATIST	VARIABLE VA FOR PARAMET r Normality geneity of es if neith iro-Wilks value .511  ************* ICS StdDev 26.56 30.98 18.79 33.69 27.63 Min	RI3 (F N RIC ANALYS of Residu variance(a er test re Levenes Test Stat 0.324  ********  StdErr 13.28 15.49 9.40 16.84 13.82  Max 91.27 118.13	>highes  A D  /TG re-analy  SIS  pals alphabsolute reserved, other  Levenes P-value 0.858  *********  Coef of V 44.42 41.40 33.76 57.39 29.88  Sof Control 125.16	st dose  ysis (ng na-level siduals) nerwise  Concli USE P. *******  Var 9:	(no sign  /mL) )  =0.01 alph non-para  usion  ARAMETRI  *******  5% Conf. 17.53, 25.54, 25.77, 5.09, 48.51, %Reduct -25.	na-level= ametric  IC TESTS  *******  Interval 102.06 124.15 85.57 112.31 136.45  ction(mea	ence

Page 80 of 86

#### EPA MRID Number 48317001

Dose4 83.75 70.47 131.96 154.66 -54.66\* PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests Analysis of Variance (ANOVA) - overall F-test Numerator df Denominator df F-stat P-value 15 1.21 Dunnett - testing each trt mean signif. different than control Williams - test assumes dose-response relationship, testing negative trend Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC Level Mean Dunnett Isotonic Williams Tukey p-values p-value mean p-value Dose1 Dose2 Dose3 Dose4 Dose5 

 Ctrl
 59.80
 .
 68.30
 .
 .

 Dosel
 74.85
 0.860
 68.30
 0.750
 .

 Dose2
 55.67
 0.999
 68.30
 0.784
 0.865

 Dose3
 58.70
 1.000
 68.30
 0.801
 0.922

 Dose4
 92.48
 0.322
 68.30
 0.812
 0.896

 1.000 0.896 0.378 0.459 \* NON-PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests Kruskal-Wallis test - equality among treatment groups Degrees of Freedom TestStat P-value 3.89 0.422 MannWhit - testing each trt median signif. different from control Jonckheere - test assumes dose-response relationship, testing negative trend Median MannWhit p-value Jonckheere p-value Level 60.45 Ctrl 67.38 0.718 Dose1 0.678 Dose2 59.01 1.000 0.442 1.000 0.426 Dose3 60.78 Dose4 83.75 0.235 0.871 DECREASING TREND TEST SUMMARY LOWEST CONCENTRATION SIGNIF. LESS THAN CONTROL Williams >highest dose (no sign. differences) Jonckheere >highest dose (no sign. differences) \* PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests Analysis of Variance (ANOVA) - overall F-test Numerator df Denominator df F-stat P-value 15 1.21 Dunnett - testing each trt mean signif. different than control Williams - test assumes dose-response relationship, testing INCREASING trend Tukey - two-sided tests, all possible comparisons, not used for NOEC or LOEC Level Mean Dunnett Isotonic Williams Tukey p-values p-value mean p-value Dose1 Dose2 Dose3 Dose4

Page 81 of 86

```
      Ctrl
      -59.80
      .
      -59.80
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .
      .

 Dose4 -92.48 0.322
                                   -92.48
                                                  0.078
                                                               0.896
                                                                         0.378 0.459
**********************
NON-PARAMETRIC ANALYSES
                                    - use alpha-level=0.05 for all tests
     Kruskal-Wallis test - equality among treatment groups
      Degrees of Freedom TestStat P-value
                                      3.89
                                                     0.422
MannWhit - testing each trt median signif. different from control
Jonckheere - test assumes dose-response relationship, testing INCREASING trend
  Level Median
Ctrl -60.45
Dosel -67.38
Dose2 -59.01
                                MannWhit p-value Jonckheere p-value
 Level
                                                  0.678
                                                                                   0.282
                                                1.000
                                                                                   0.558
  Dose3
              -60.78
                                                1.000
                                                                                   0.574
                                                0.235
  Dose4 -83.75
                                                                                   0.129
 INCREASING TREND TEST SUMMARY LOWEST CONCENTRATION SIGNIF. GREATER THAN
CONTROL
    Williams
                                                                >highest dose (no sign. differences)
    Jonckheere
                                                                >highest dose (no sign. differences)
test for fish screen study - TEST DATA 2 4 D
ANALYSIS RESULTS FOR VARIABLE VAR14 ( M testosterone (ng/mL) )
TESTS OF ASSUMPTIONS FOR PARAMETRIC ANALYSIS
Shapiro-Wilks test for Normality of Residuals -- alpha-level=0.01
Levenes test for homogeneity of variance(absolute residuals) -- alpha-level=0.05
Use parametric analyses if neither test rejected, otherwise non-parametric
analyses.
  Shapiro-Wilks Shapiro-Wilks Levenes Levenes Conclusion
     Test Stat P-value Test Stat P-value
                                                              0.858 USE PARAMETRIC TESTS
                                               0.324
        0.958
                          0.511
*************************
BASIC SUMMARY STATISTICS
 Level N Mean StdDev StdErr Coef of Var 95% Conf.Interval
                           •
  Ctrl 0
                                                .
  Dosel 0
  Dose2 0
                                               .
  Dose3 0
                                                .
                                 .
  Dose4 0
            Median Min Max %of Control(means) %Reduction(means)
 Level
  Ctrl
  Dose1
  Dose2
```

Page 82 of 86

```
Dose3
 Dose4
**************************
PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
   Analysis of Variance (ANOVA) - overall F-test
   Numerator df Denominator df F-stat P-value
                 15
                                1.21
                                          0.346
MannWhit - testing each trt median signif. different from control
Jonckheere - test assumes dose-response relationship, testing negative trend
       Median MannWhit p-value Jonckheere p-value
Level
 Ctrl
 Dose1
 Dose2
 Dose3
 Dose4
  Jonckheere
*************************
PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
   Analysis of Variance (ANOVA) - overall F-test
    Numerator df Denominator df F-stat
                                          P-value
                                          0.346
                 15
                                1.21
MannWhit - testing each trt median signif. different from control
Jonckheere - test assumes dose-response relationship, testing INCREASING trend
       Median MannWhit p-value Jonckheere p-value
Level
 Ctrl
 Dose1
 Dose2
 Dose3
 Dose4
  Jonckheere
                                      Dose1
test for fish screen study - TEST DATA 2 4 D
ANALYSIS RESULTS FOR VARIABLE VAR15 ( F 17b-estradiol (ng/mL) )
TESTS OF ASSUMPTIONS FOR PARAMETRIC ANALYSIS
Shapiro-Wilks test for Normality of Residuals -- alpha-level=0.01
Levenes test for homogeneity of variance(absolute residuals) -- alpha-level=0.05
Use parametric analyses if neither test rejected, otherwise non-parametric
analyses.
 Shapiro-Wilks Shapiro-Wilks Levenes
                                     Levenes Conclusion
   Test Stat P-value Test Stat P-value
               0.511
                          0.324 0.858 USE PARAMETRIC TESTS
************************
BASIC SUMMARY STATISTICS
Level N Mean StdDev StdErr Coef of Var 95% Conf.Interval Ctrl 0 . . . . . . . . . . . . . . . .
```

Page 83 of 86

```
Dosel 0
 Dose2 0
 Dose3 0
 Dose4 0
Level Median Min Max %of Control(means) %Reduction(means)
 Ctrl
 Dose1
 Dose2
 Dose3
 Dose4
*******************
PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
   Analysis of Variance (ANOVA) - overall F-test
   Numerator df Denominator df F-stat P-value 4 15 1.21 0.346
MannWhit - testing each trt median signif. different from control
Jonckheere - test assumes dose-response relationship, testing negative trend
Level Median
                     MannWhit p-value
                                           Jonckheere p-value
 Ctrl
 Dose1
 Dose2
 Dose3
 Dose4
  Jonckheere
                                      Dose1
************************
PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests
   Analysis of Variance (ANOVA) - overall F-test
    Numerator df Denominator df F-stat P-value
                 15
                                1.21
                                          0.346
MannWhit - testing each trt median signif. different from control
Jonckheere - test assumes dose-response relationship, testing INCREASING trend
                    MannWhit p-value Jonckheere p-value
Level
       Median
 Ctrl
 Dose1
 Dose2
 Dose3
 Dose4
  Jonckheere
                                      Dose1
test for fish screen study - TEST DATA 2 4 D
ANALYSIS RESULTS FOR VARIABLE VAR16 ( M 17b-estradiol (ng/mL) )
TESTS OF ASSUMPTIONS FOR PARAMETRIC ANALYSIS
Shapiro-Wilks test for Normality of Residuals -- alpha-level=0.01
Levenes test for homogeneity of variance(absolute residuals) -- alpha-level=0.05
```

Use parametric analyses if neither test rejected, otherwise non-parametric Shapiro-Wilks Shapiro-Wilks Levenes Levenes Conclusion Test Stat P-value Test Stat P-value 0.958 0.511 0.324 0.858 USE PARAMETRIC TESTS \* BASIC SUMMARY STATISTICS Level N Mean StdDev StdErr Coef of Var 95% Conf.Interval Ctrl 0 Dose1 0
Dose2 0
Dose3 0 Dose4 0 Level Median Min Max %of Control(means) %Reduction(means) Ctrl . . Dose1 Dose2 Dose3 Dose4 \* PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests Analysis of Variance (ANOVA) - overall F-test P-value Numerator df Denominator df F-stat 15 1.21 MannWhit - testing each trt median signif. different from control Jonckheere - test assumes dose-response relationship, testing negative trend Level Median MannWhit p-value Jonckheere p-value Ctrl Dose1 Dose2 Dose3 Dose4 Jonckheere Dose1 \* PARAMETRIC ANALYSES - use alpha-level=0.05 for all tests Analysis of Variance (ANOVA) - overall F-test Numerator df Denominator df F-stat P-value 1.21 MannWhit - testing each trt median signif. different from control Jonckheere - test assumes dose-response relationship, testing INCREASING trend MannWhit p-value Jonckheere p-value Level Median Ctrl Dose1 Dose2 Dose3

Page 85 of 86

# Data Evaluation Record on the Fish Short-Term Reproduction Assay with 2,4-D

Dose4 . . . . . . . . . . Jonckheere Dose1

# DATA EVALUATION RECORD

2,4-DICHLOROPHENOXY ACETIC ACID (2,4-D)

Study Type: OCSPP 890.1550, Steroidogenesis Assay

EPA Contract No. EP10H001452 Task Assignment No. 2-26-2012 (MRID 48614305)

Prepared for
Health Effects Division
Office of Pesticide Programs
U.S. Environmental Protection Agency
2777 South Crystal Drive
Arlington, VA 22202

Prepared by CSS-Dynamac Corporation 1910 Sedwick Road, Building 100, Suite B Durham, NC 27713

Primary Reviewer:	Signature:	Konnie J. Bever G.
Ronnie J. Bever Jr., Ph.D., D.A.B.T.	Date:	3/14/2012
Secondary Reviewer	Signature:	Lutenberg
Scott D. Studenberg, Ph.D., D.A.B.T.	Date:	3/22/2012
Program Manager:	Signature:	Jack D. Eny
Jack D. Early, M.S.	Date:	3/27/2012
Quality Assurance:	Signature:	Jack D. Eny
Jack D. Early, M.S.	Date:	3/27/2012

This Data Evaluation Record may have been altered by the Health Effects Division subsequent to signing by CSS-Dynamac Corporation personnel.

The US EPA Endocrine Disruptor Screening Program (EDSP) Tier 1 screening battery is comprised of eleven screening assays intended to identify a chemical's likely endocrine bioactivity, i.e., its potential to interact with the estrogen, androgen, or thyroid (E, A, or T) pathways. The robustness of the Tier 1 battery is based on the strengths of each individual assay to identify potential endocrine bioactivity with complementary endpoints within the assay, where available, and redundancy across the battery. Thus, the results of each individual assay should not be considered in isolation but rather should be considered in the context of other assays in the battery as well as Other Scientifically Relevant Information (OSRI). In order to determine if a chemical has the potential to interact with the E, A or T pathways, a Weight of Evidence (WoE) evaluation of Tier 1 assay results, in combination with the findings in the OSRI, should be undertaken (refer to the WoE Document).

Primary Reviewer: Gregory Akerman

**Health Effects Division** 

Secondary Reviewer: John Liccione, Ph.D.

**Health Effects Division** 

Signature:

Signature:

Date

Template version 08/2011

# DATA EVALUATION RECORD

STUDY TYPE: Steroidogenesis Assay (H295R Cells); OCSPP 890.1550

PC CODE: 030001 DP BARCODE: D398637

TXR#: 0052104 CAS No.: 94-75-7

TEST MATERIAL (PURITY): 2,4-D (98.5% a.i.)

SYNONYMS: 2,4-dichlorophenoxyacetic acid; 2,4-D acid, 2-(2,4-dichlorophenoxy)acetic

acid

**CITATION:** LeBaron, M.J.; Kan, H.L; Perala, A.W. (2011) Evaluation of 2,4-

dichlorophenoxyacetic acid (2,4-D) in the *in vitro* steroidogenesis assay. Toxicology & Environmental Research and Consulting, The Dow Chemical Company, Midland, MI. Laboratory Study No.: 111038, October 20, 2011.

MRID 48614305. Unpublished.

**SPONSOR:** Industry Task Force II on 2,4-D Research Data, 1900 K Street NW,

Washington, DC

TEST ORDER #: CON-030001-1

**EXECUTIVE SUMMARY:** In a steroidogenesis assay (MRID 48614305), H295R cells cultured *in vitro* in 24-well plates were incubated with 2,4-D (98.5% purity, Lot # 2006 2433 8006-USA) at log concentrations of 10<sup>-10</sup> to 10<sup>-4</sup> M for 48 hours in triplicate in three independent experiments. Dimethyl sulfoxide (DMSO) was used as the vehicle, at a final concentration in the assay of 0.1%.

Testosterone and estradiol levels were measured using LC-APPI-MS/MS. A Quality Control (QC) plate was run concurrently with each independent run of a test chemical plate to demonstrate that the assay responded properly to positive control agents at two concentration levels. The positive controls included a known inhibitor (prochloraz) and a known inducer (forskolin) of estradiol and testosterone production.

Guideline acceptability recommendations and requirements were generally met, including lack of cytotoxicity, adequate production of testosterone and estradiol, acceptable reproducibility (low %CV), and appropriate induction and inhibition with positive controls, with two exceptions. The required concentration for estradiol production (40 pg/mL) was not met in any run for blank or solvent control, with values ranging from 26.7 to 38.3 pg/mL, but the basal concentrations were greater than 2.5-fold of the minimum detection level (10 pg/mL) as required. Secondly, all three

runs of 1  $\mu$ M prochloraz only reduced estradiol to 0.6-fold that of the solvent control (instead of 0.5-fold).

2,4-D had no effect on testosterone production at concentrations up to  $10^{-4}$  M, and no effect on estradiol production at concentrations up to  $10^{-5}$  M. At the highest concentration tested ( $10^{-4}$  M), 2,4-D increased (p≤0.05) estradiol levels by 20% (1.2-fold) in all three runs relative to the DMSO-treated cells.

Based on hormone responses in each of the independent runs, 2,4-D treatment resulted in statistically significant and reproducible increases in estradiol production. 2,4-D treatment did not result in statistically significant and reproducible alterations in testosterone production.

The assay **satisfies** the EDSP Tier 1 Test Order requirements for a Steroidogenesis Assay (OCSPP 890.1550).

**<u>COMPLIANCE</u>**: Signed and dated GLP Compliance, Quality Assurance, and Data Confidentiality statements were provided.

#### I. MATERIALS AND METHODS

### **MATERIALS**

Toxicology & Environmental Research and Consulting, **Test Facility:** 

The Dow Chemical Company

**Location:** Midland, MI **Study Director:** H. L. Kan

Other Personnel: M. J. LeBaron (lead scientist), A. W. Perala (analytical chemist), B. B.

Gollapudi.(technical reviewer)

March 14, 2011 - October 20, 2011 **Study Period:** 

2,4-D **Test Substance:** 

> **Description:** Off-white powder

2006 2433 8006-USA (not provided) Lot # (expiration date):

**Purity:** 98.5%

Soluble in DMSO up to 0.1 M; up to 315 mg/L in water **Solubility:** 

Vapor pressure:  $1.9 \times 10^{-5}$  Pa at  $25^{\circ}$ C

**Stability:** 2 year shelf life; stability under test conditions not provided

**Storage conditions:** Ambient 94-75-7 CAS #: 221.0 Molecular weight:

Structure:

Forskolin **Positive Control (Inducer):** 

> **Description (molecular weight):** White powder (410.50)

Sigma-Aldrich (St. Louis, MO) Source: Lot # (expiration date): 097K50653 (not provided)

**Purity:** 99%

**Solubility (in solvent):** Soluble in DMSO up to 0.01 M

Storage conditions: Ambient CAS#: 66575-29-9

**Positive Control (Inhibitor):** Prochloraz

> **Description (molecular weight):** Off-white powder (376.67) Source: Sigma-Aldrich (St. Louis, MO)

Lot # (expiration date): SZE6220X (not provided)

**Purity:** 99.1%

Soluble in DMSO up to 0.01 M **Solubility (in solvent):** 

Storage conditions: Ambient 67747-09-5 CAS #:

5. Solvent/Vehicle Control: Dimethyl Sulfoxide (DMSO)

**Description (molecular weight):** Clear liquid (78.13)

Sigma-Aldrich (St. Louis, MO) Source:

Lot # (expiration date): 16496CPV and 68996LMV (not provided)

≥99.9% **Purity:** 

**Storage conditions:** Room temperature

67-68-5 CAS#:

Not provided. DMSO is a Guideline recommended solvent in which 2,4-D is Justification for choice of solvent:

reasonably soluble.

0.1% v/vFinal concentration in assay:

**6. Stock Medium:** Dulbecco's modified Eagle's medium/Ham's F12

nutrient mixture

Source: Sigma-Aldrich (St. Louis, MO)
Lot # (expiration date): RNBB 6106 (not provided)
Sodium bicarbonate: Component of stock medium

Nu-Serum: 2.5%, BD Biosciences; Lot # 81515; tested for background hormone

concentrations by performing laboratory BD Biosciences; Lot # 88964 & 05245

ITS+premix: BD Biosciences; Lot # 88964 & 05245
Other components: 2.5 mM L-glutamine (Life Technologies, GIBCO, Grand Island, NY), and

penicillin/streptomycin (25 IU/mL/25 µg/mL) (GIBCO).

7. <u>Test Cells</u>: H295R human adrenocortical carcinoma cells (ATCC CLR-2128; Batch # not provided) at passage 7.5 – 8.5 were incubated in the stock medium. Incubation conditions were at approximately 37° C with 5% CO<sub>2</sub> for approximately 24 hours prior to exposure.

The following performance criteria were met (indicated by an "x"):

X	Cell passage identifier. Cell Passage #: 7.5 – 8.5
X	Cells frozen down at passage 5
X	Frozen cells cultured for 4 additional passages
X	Total number of passages does not exceed 10

# **B.** METHODS

# 1. Pre-Test Information

- **a.** <u>Hormone Assay Interference Test</u>: A separate chemical interference assay was not performed, as LC/APPI-MS/MS was used along with internal spike-in controls to measure extraction and quantitation efficiency.
- **b.** Hormone Extraction: Briefly, standards and samples were vortexed with methylene chloride containing both testosterone and estradiol internal standards, and the organic phase were transferred to clean vials. The aliquots were evaporated to dryness, and the final residue in each vial was derivatized with dansyl chloride in sodium bicarbonate buffer (100mM, pH 10.5). Samples were then subjected to LC/APPI-MS/MS analysis.
- **c.** <u>Laboratory Proficiency Test</u>: The Sponsor stated that Laboratory proficiency assays were performed to optimize and validate the H295R steroidogenesis assay. These non-GLP, unpublished experimental results with positive control chemicals demonstrated laboratory proficiency and validation and are reported separately (data not presented and report not cited).
- 2. <u>Test Solutions</u>: 2,4-D, forskolin, and prochloraz were dissolved in DMSO and diluted 1:1000 in the final treatment medium. When added to the cell culture plates, these dilutions yielded final concentrations of 1 μM and 10 μM for forskolin, 0.1 μM and 1 μM for prochloraz, and 10<sup>-4</sup>, 10<sup>-5</sup>, 10<sup>-6</sup>, 10<sup>-7</sup>, 10<sup>-8</sup>, 10<sup>-9</sup>, and 10<sup>-10</sup> M for 2,4-D. The selected concentrations of the test material in the dosing solutions used for the treatment in a definitive assay (generally the first assay) were verified analytically. The final concentration of DMSO in the medium was held constant at 0.1% (v/v).

- 3. Cell Plating and Preincubation: Cells were maintained in the stock medium described above. The hormone concentrations in undiluted Nu-Serum were 974 pg/mL for testosterone and 681 pg/mL for estradiol. H295R cells were grown for five passages, frozen in liquid nitrogen, then thawed and cultured for seven to eight additional passages prior to use in the assay. The cells were plated into wells of a 24-well cell culture plate at a density of approximately 200,000 300,000 cells in 1 mL (50 60% confluency). The cells were then placed into an incubator (approximately 37° C/5% CO<sub>2</sub>) for approximately 24 hours prior to chemical exposure. After a 24-hour pre-incubation, the cells were checked microscopically for good attachment and proper morphology.
- **Exposure:** The medium was removed from the cells and replaced with 1 mL of medium containing concentrations of 2,4-D (or solvent) in triplicate, according to the schematic presented in Table 1.

TABLE 1.	Dosing Schen	Dosing Schematic for the Exposure of H295R Cells to 2,4-D (Final Concentrations in M). <sup>a</sup>									
	1	1 2 3 4 5 6									
A	DMSO	DMSO	DMSO	$10^{-7}$	$10^{-7}$	$10^{-7}$					
В	$10^{-4}$	$10^{-4}$	$10^{-4}$	$10^{-8}$	$10^{-8}$	$10^{-8}$					
С	$10^{-5}$	$10^{-5}$	$10^{-5}$	$10^{-9}$	$10^{-9}$	$10^{-9}$					
D	$10^{-6}$	$10^{-6}$	$10^{-6}$	$10^{-10}$	$10^{-10}$	$10^{-10}$					

a Data were obtained from page 16 of the study report.

A concurrent QC plate was included with each of the three independent runs of the test chemical plates to demonstrate the assay's response to forskolin (an inducer of testosterone and estradiol production) and prochloraz (an inhibitor of testosterone and estradiol production). The QC plate was prepared and dosed in the same manner with either forskolin or prochloraz, according to the schematic presented in Table 2.

TA	TABLE 2. Dosing Schematic for the QC Plate for Positive Controls (Final Concentrations in μM). <sup>a</sup>									
	1 2		3	4	5	6				
A	Blank <sup>b</sup>	Blank	Blank	Blank + methanol <sup>c</sup>	Blank + methanol	Blank + methanol				
B DMSO		DMSO	DMSO	DMSO + methanol	DMSO + methanol	DMSO + methanol				
C	Forskolin	Forskolin	Forskolin	Prochloraz	Prochloraz	Prochloraz				
1 μΜ 1 μΜ		1 μM	1 μM	0.1 μΜ	0.1 μΜ	0.1 μΜ				
<b>D</b> Forskolin		Forskolin	Forskolin	Prochloraz	Prochloraz	Prochloraz				
	10 μΜ	10 μΜ	10 μΜ	1 μM	1 μM	1 μΜ				

- a Data were obtained from page 15 of the study report.
- b Blank wells received medium only.
- c Methanol (MeOH) was added after the exposure was terminated and the medium was removed.

Following dosing, the plates were incubated for approximately 48 hours under the conditions previously described. After the incubation period, two aliquots of medium were collected and frozen at  $-80^{\circ}$ C until further processing.

5. <u>Cell Viability/Cytotoxicity Assay</u>: The cell viability/cytotoxicity testing was conducted in the QC plate and in the test chemical exposure plate immediately after termination of the exposure experiments. In addition to viability/cytotoxicity testing, cells were checked for the degree of confluence, homogeneity from well-to-well, and any signs of cytotoxicity or

altered morphology. A subjective parallel visual assessment of each well was conducted. The CellTiter 96 Aqueous One Solution Cell Proliferation Assay kit (Madison Wisconsin, Cat# G3580) was used for assessment of cell viability testing, as per manufacturer's instructions. This cell proliferation assay kit is a colorimetric modified MTT cell viability assay. The Aqueous One Solution contains a tetrazolium compound [3-(4,5-dimethylthiazol-2-yl)-5(3-carboxymethoxyphenyl)-2-(4-sulfophenyl)-2H-tetrazolium, inner salt; MTS] and an electron coupling reagent (phenazine ethosulfate; PES). The quantity of formazan product as measured by the amount of 490 nm absorbance is directly proportional to the number of viable cells in culture. The minimum cell viability required per well was ≥80%; wells with lower viability were excluded in the final data analysis.

**6.** Hormone Measurement System: The concentration of each hormone (testosterone and estradiol) was evaluated in triplicate in the culture medium. Hormones were extracted from the samples with methylene chloride, derivatized with dansyl chloride, and analyzed by LC/APPI-MS/MS. The lower limit of quantification (LLQ) for this measurement system was 25 pg/mL for testosterone and 10 pg/mL for estradiol.

The following performance criteria were met (indicated by an "x"):

X	Method detection limit (100 pg/mL testosterone; 10 pg/mL estradiol)
X	Spiked sample recovery acceptable for two concentrations of testosterone and estradiol (mean measured amount
Λ	from triplicate samples within 30% of nominal concentration)
NA	Hormone cross-reactivity (antibody-based assays only; ≤30% of basal production of the respective hormone)
X	Solvent control within 75% range below maximum response on standard curve
NA	Test compound tested for interference with measurement system

C. <u>DATA ANALYSIS</u>: To evaluate the relative increase or decrease of hormone production after test chemical exposure, the results were normalized to the mean solvent control (SC) value for each assay and results were expressed as fold-change relative to the SC in each exposure plate. All data were expressed as mean ± Standard Deviation (SD). Relative changes were calculated as follows:

Relative Change = Hormone concentration in each well ÷ Hormone concentration of mean solvent (vehicle) control

Prior to conducting statistical analyses, the assumptions of normality and variance homogeneity were evaluated. Homogeneity of variance were evaluated by Bartlett's test and normality by Shapiro-Wilk's test at alpha = 0.01. If the data were not homogeneous or normally distributed, then the data were transformed to approximate homogeneity or a normal distribution. If the data were homogeneous and approximately normally distributed, differences between chemical treatments and SC were analyzed using a parametric analysis of variance followed by Dunnett's test, if significant. If the data were not homogeneous or normally distributed, a non-parametric test was used (Kruskal-Wallis) and if significant, was followed by the Wilcoxon rank sum test with a Bonferroni-Holm correction. Differences were considered significant at  $p \le 0.05$ . These statistical analyses were considered appropriate.


#### II. RESULTS


- **A.** TEST COMPOUND: Precipitation of the test compound was not reported at any concentration. Cytotoxicity was not observed at any tested 2,4-D concentration. The variability (CV) between the runs for the solvent controls was 4.6% (testosterone) and 14.6% (estradiol), both within the guideline recommendation of 30%. The variability (CV) within each run for the solvent controls in the testosterone and estradiol assays was 13.6% or less, which was within the guideline recommendation of 30%.
  - 2,4-D had no effect on testosterone levels at any concentration tested (Table 3, Figure 1), and a slight effect on estradiol levels at 2,4-D concentrations of 10  $\mu$ M or less. However, at a concentration of 100  $\mu$ M, 2,4-D increased (p $\leq$ 0.05) estradiol levels by 20% relative to DMSO-treated cells in all three runs (Figure 2).

TABLE 3. Mean (±SD) Hormone Concentrations Following Treatment with 2,4-D for 48 Hours. <sup>a</sup>										
Nominal Concentration (µM)	Run 1	Run 2	Run 3	Run 1	Run 2	Run 3	Mean	± SD	Statistical Significance <sup>b</sup>	
	Testo	sterone (p	g/mL)		Fold Difference					
DMSO	842.0	735.0	778.0	_	_	_	_	_		
0.0001	776.0	691.7	746.3	0.9	0.9	1.0	0.9	0.0	_	
0.001	805.7	696.0	752.0	1.0	0.9	1.0	1.0	0.0	_	
0.01	737.7	696.3	726.7	0.9	0.9	0.9	0.9	0.0		
0.1	752.0	647.0	745.3	0.9	0.9	1.0	0.9	0.0	_	
1	796.7	665.0	734.0	0.9	0.9	0.9	0.9	0.0	_	
10	787.7	684.3	747.0	0.9	0.9	1.0	0.9	0.0	_	
100	714.3	610.3	728.0	0.8	0.8	0.9	0.9	0.1	_	
	Est	radiol (pg/	/mL)		-	Fold	Difference		-	
DMSO	29.8	35.7	37.0		_	_	_	_	_	
0.0001	24.7	37.2	36.0	0.8	1.0	1.0	0.9	0.1	_	
0.001	24.7	35.7	34.2	0.8	1.0	0.9	0.9	0.1		
0.01	25.6	34.8	37.8	0.9	1.0	1.0	1.0	0.1	_	
0.1	26.5	36.4	34.9	0.9	1.0	0.9	1.0	0.1	_	
1	25.3	33.7	35.1	0.8	0.9	0.9	0.9	0.1		
10	27.0	33.8	37.9	0.9	0.9	1.0	1.0	0.1	_	
100	34.6	44.3	45.1	1.2	1.2	1.2	1.2	0.0	Runs 1, 2, 3	

a Data were obtained from page 28 of the study report. The lower limit of quantification (LLQ) was 25 pg/mL for testosterone and 10 pg/mL for estradiol.

FIGURE 1. Change in Testosterone Production Relative to 2,4-D Concentration.

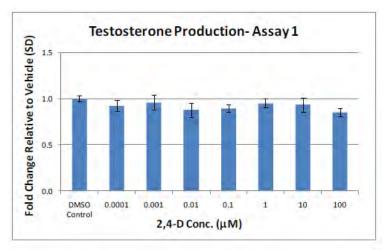

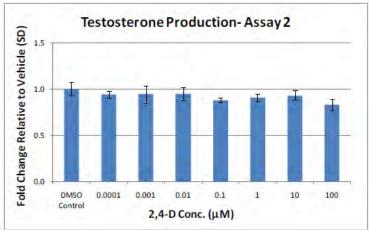




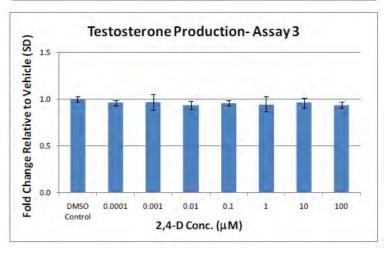


FIGURE 2. Change in Estradiol Production Relative to 2,4-D Concentration.

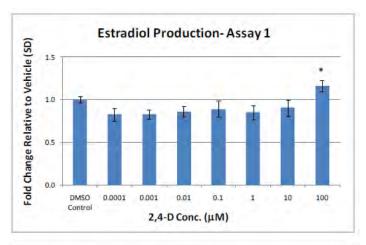

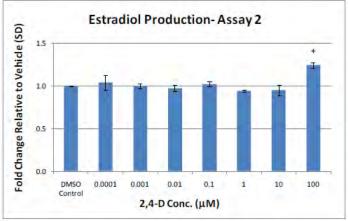

\* Estradiol production significantly different from solvent control (p<0.05)

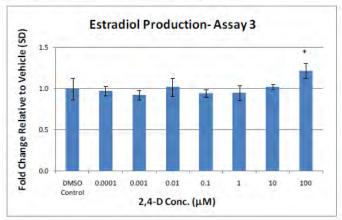

<sup>\*</sup> Estradiol production significantly different from solvent control (p<0.05)

**B.** <u>CYTOTOXICITY</u>: 2,4-D was not cytotoxic to the H295R cells (Table 4). Viability in each run ranged from 92.5–133% relative to solvent controls.

TABLE 4. Mean (±SD) MTT Cell Viability Results after Treatment with Forskolin, Prochloraz, or 2,4-D for 48 Hours. <sup>a</sup>							
,	Concentration		<b>Tiability</b>				
Compound	(μΜ)	Mean	SD				
Media	NA	96.6	1.59				
Media + Methanol	NA	29.4	2.35				
DMSO + Methanol	NA	30.9	3.18				
Forskolin	1	121.7	9.91				
Forskolin	10	114.1	5.25				
Prochloraz	0.1	110.1	9.50				
Prochloraz	1	99.9	7.40				
2,4-D	0.0001	116.6	8.05				
2,4-D	0.001	110.9	5.37				
2,4-D	0.01	111.7	8.35				
2,4-D	0.1	104.9	4.92				
2,4-D	1	106.4	5.89				
2,4-D	10	110.7	5.16				
2,4-D	100	106.9	5.12				

a Values were calculated by the reviewers from data obtained from page 41 of the study report.

C. QC PLATE: The hormone concentrations after exposure to the reference chemicals, SC, and blank samples, as well as the fold difference change relative to SC (individual and mean ± SD) for the three assay runs are presented in Table 5. The minimum basal hormone production levels (500 pg/mL for testosterone, 40 pg/mL for estradiol) were met in both blank and SC wells for testosterone (Table 6). The required concentration for estradiol was not met in any run for blank or SC, with values ranging from 26.7–38.3 pg/mL. Forskolin at 10 μM induced testosterone an average of 2.4-fold and estradiol 20-fold over the SC. Prochloraz at 1 μM inhibited concentrations of testosterone by 0.3-fold and estradiol by 0.6-fold compared to SC. Guideline requirements for basal hormone induction and inhibition were met for testosterone and estradiol, except that all three runs of 1 μM prochloraz only reduced estradiol to 0.6-fold that of the SC rather than 0.5-fold.

The variability (%CV) between the runs (calculated by the reviewer) based on the absolute hormone concentrations in the SC were 4.6% for testosterone and 14.6% for estradiol, and within the guideline criteria of  $\leq$ 30% for the assays. The %CVs within each run for the QC plates were 2.5-6.6% for testosterone and 3.8-10.1% for estradiol, which were within the guideline criteria ( $\leq$ 30%).

TABLE 5. Mean (±SD) Hormone Concentrations Following Treatment with Forskolin or Prochloraz for 48 Hours. <sup>a</sup>								
Concentration	Run 1	Run 2	Run 3	Run 1	Run 2	Run 3	Mean	± SD
(μΜ)	Testosterone (pg/mL)			Fold Difference (Relative to DMSO)				
Blank	1041	936	1040		_			_
DMSO	778	719	781					
1 μM Forskolin	1377	1217	1310	1.8	17	1.7	1.7	0.0
10 μM Forskolin	2007	1757	1733	2.6	2.4	2.2	2.4	0.2
0.1 µM Prochloraz	515	514	534	0.7	0.7	0.7	0.7	0.0
1 μM Prochloraz	250	233	218	0.3	0.3	0.3	0.3	0.0
	Est	tradiol (pg/n	ıL)	Fold Difference (Relative to DMSO)				
Blank	26.3	34.8	37.2	_	_	_	_	_
DMSO	28.5	34.2	38.3	_	_			1
1 μM Forskolin	348	366	492	12.2	10.7	12.8	11.9	1.1
10 μM Forskolin	610	626	780	21.4	18.3	20.4	20.0	1.6
0.1 μM Prochloraz	26.4	29.7	28.6	0.9	0.9	0.5	0.8	0.3
1 μM Prochloraz	17.6	19.5	21.9	0.6	0.6	0.6	0.6	0.0

Data were obtained from page 29 of the study report. The lower limit of quantification (LLQ) was 25 pg/mL for testosterone and 10 pg/mL for estradiol.

#### III. DISCUSSION AND CONCLUSIONS

- **A.** <u>INVESTIGATOR'S CONCLUSIONS</u>: Based on the combined hormone responses in each of three independent H295R steroidogenesis assays, it was determined that 2,4-D treatment resulted in a statistically significant increase in estradiol production at the assay limit-concentration of 100 μM (10<sup>-4</sup> M), while no effects on estradiol production were seen at lower concentrations. There were no effects on testosterone production at any concentration of 2,4-D. The slight (1.2-fold or less) increase in estradiol production did not meet the 1.5-fold cut-off criteria established in the validation program for this assay, and hence was not interpreted to be biologically relevant.
- **B.** AGENCY COMMENTS: Guideline acceptability recommendations and requirements were generally met, including lack of cytotoxicity, adequate production of testosterone and estradiol, acceptable reproducibility (low %CV), and appropriate induction and inhibition with positive controls. The exceptions included: the required concentration for estradiol production (40 pg/mL) was not met in any run for blank or solvent control, with values ranging from 26.7–38.3 pg/mL; and all three runs of 1 μM prochloraz only reduced estradiol to 0.6-fold that of the solvent control rather than 0.5-fold.
  - 2,4-D had no effect on testosterone production at concentrations up to  $10^{-4}$  M, and no effect on estradiol production at concentration up to  $10^{-5}$  M. However, at the highest concentration tested ( $10^{-4}$  M), 2,4-D increased ( $p \le 0.05$ ) estradiol levels by 20% in all three runs relative to the DMSO-treated cells.

Based on hormone responses in each of the independent runs, 2,4-D treatment resulted in a statistically significant and reproducible alteration in estradiol production.

- C. <u>STUDY DEFICIENCIES</u>: The following deficiencies were noted that are not considered to have had an adverse impact on the results, interpretation or conclusions of this study:
  - Stability of 2,4-D in the test system was not reported.
  - The report erroneously specified an LLQ of 25 mg/mL. The assay must be able to quantify 500 pg/mL according to the Guideline, and this assay was able to do so.
  - There were two minor departures from the Guideline: (i) the required concentration for estradiol was not met in any run for blank or SC, ranging from 26.7 to 38.3 pg/mL rather than 40 pg/mL; and (ii) all three runs of 1 μM prochloraz only reduced estradiol to 0.6-fold that of the SC rather than 0.5-fold.